

ภาคผนวกที่ 6

เอกสารการสอบเทียบความถูกต้อง
ของเครื่องมือ

**ตารางสรุปรายการเอกสารการสอบเทียบความถูกต้องของเครื่องมือเก็บตัวอย่าง
และเครื่องมือตรวจวิเคราะห์คุณภาพสิ่งแวดล้อม**

รายการตรวจวัด	เครื่องมือเก็บตัวอย่าง	เครื่องมือตรวจวิเคราะห์
	ชื่อเครื่องมือ	ชื่อเครื่องมือ
1. คุณภาพอากาศในบรรยากาศ - TSP	- High Volume Air Sampler & Blower -No. B16, B16 -No. B19, B19 -No. B29, B29 -No. B33, B33	- Electronic Balance
- PM-10	- High Volume PM-10 Air Sampler & Blower -No. B6, B6 -No. B10, B10 -No. B18, B18 -No. B20, B20	- Electronic Balance
2. ระดับเสียงในบรรยากาศ - Leq 24 hr, Lmax	- Acoustic Calibrator - Sound Level Meter - ACO-B06 - ACO-B07 - ACO-B17 - ACO-B21 - ACO-B40	-
3.การตรวจวิเคราะห์คุณภาพน้ำ 1. pH 2. Turbidity 3. Total Suspended Solids 4. Total Dissolved Solids 5. Sulfate 6. Total Hardness 7. Total Iron 8. Lead 9. Arsenic 10.Cadmium	- - - - - - - - - -	- pH Meter - Turbidity Meter - Electronic Balance - Electronic Balance - Spectrophotometer - Electronic Balance - Inductively Coupled Plasma (ICP) - Inductively Coupled Plasma (ICP) - Hydride Generation/ Atomic Absorption Spectrophotometer - Inductively Coupled Plasma (ICP)



บริษัท เอส.พี.เอส. คอนซัลติ้ง เซอร์วิส จำกัด
S.P.S. CONSULTING SERVICE CO., LTD.
7 ซอยพหลโยธิน 24 ถนนพหลโยธิน แขวงจอมพล เขตจตุจักร กรุงเทพฯ 10900
7 Soi Phaholyothin 24, Phaholyothin Rd., Jompol, Chatuchak, Bangkok 10900
Tel : (662) 939-4370-72 Fax : (662) 513-4221 E-mail : sale@spscon.com, www.spscon.com

High Volume Air Sampler Calibration Report

Calibration Method : Multipoint Orifice Flow Transfer Standard

Model : TE 5025A

S/N : 3095

Calibration Data

High Volume Air Sampler Data		Calibration Data		
Recorder No.	Blower No.	Date	Actual Flowrate (ft ³ /min)	R ²
B01	B01	09/02/2022	y = 1.255x-7.443	0.998
B02	B02	02/02/2022	y = 1.075x+1.871	0.999
B03	B03	04/02/2022	y = 1.032x+1.126	0.997
B04	B04	04/02/2022	y = 1.158x-3.770	0.995
B05	B05	02/02/2022	y = 1.199x-5.374	1.000
B06	B06	01/02/2022	y = 1.215x-6.623	0.995
B07	B07	01/02/2022	y = 1.142x-4.465	0.997
B08	B08	02/02/2022	y = 1.241x-8.074	0.999
B09	B09	08/02/2022	y = 1.206x-5.652	0.995
B10	B10	07/02/2022	y = 1.095x+0.184	0.998
B11	B11	10/02/2022	y = 1.099x-2.021	0.996
B12	B12	09/02/2022	y = 1.169x-3.784	1.000
B13	B13	03/02/2022	y = 1.163x-4.662	0.996
B14	B14	07/02/2022	y = 1.169x-3.363	0.998
B15	B15	03/02/2022	y = 1.106x-1.273	0.998
B16	B16	09/02/2022	y = 1.218x-6.757	0.997
B17	B17	07/02/2022	y = 1.132x-1.890	0.998
B18	B18	16/02/2022	y = 1.239x-7.560	0.999
B19	B19	16/02/2022	y = 1.265x-8.934	0.997
B20	B20	03/02/2022	y = 1.199x-6.304	0.998
B21	B21	17/02/2022	y = 1.120x-2.616	0.997
B22	B22	08/02/2022	y = 1.216x-6.597	0.995
B23	B23	03/02/2022	y = 1.139x-3.341	0.999
B24	B24	03/02/2022	y = 1.126x-2.172	1.000
B25	B25	09/02/2022	y = 1.016x+2.185	0.996
B26	B26	04/02/2022	y = 1.122x-2.540	0.997
B27	B27	08/02/2022	y = 1.192x-6.584	0.997
B28	B28	04/02/2022	y = 1.254x-8.360	0.995
B29	B29	02/02/2022	y = 1.217x-6.791	0.996
B30	B30	04/02/2022	y = 1.162x-4.303	0.997
B31	B31	16/02/2022	y = 1.101x-0.556	0.998
B32	B32	04/02/2022	y = 1.208x-5.034	0.997
B33	B33	07/02/2022	y = 1.242x-5.616	0.999
B34	B34	09/02/2022	y = 1.240x-8.273	0.999

Calibrated by :



บริษัท เอส.พี.เอส. คอนซัลติ้ง เซอร์วิส จำกัด
S.P.S. CONSULTING SERVICE CO., LTD.
7 ซอยพหลโยธิน 24 ถนนพหลโยธิน แขวงจอมพล เขตจตุจักร กรุงเทพฯ 10900
7 Soi Phaholyothin 24, Phaholyothin Rd., Jompol, Chatuchak, Bangkok 10900
Tel : (662) 939-4370-72 Fax : (662) 513-4221 E-mail : sale@spscon.com, www.spscon.com

High Volume Air Sampler Calibration Report

Calibration Method : Multipoint Orifice Flow Transfer Standard


Model : TE 5025A

S/N : 3095

Calibration Data

High Volume Air Sampler Data		Calibration Data		
Recorder No.	Blower No.	Date	Actual Flowrate (ft ³ /min)	R ²
B35	B35	16/02/2022	y = 1.274x-9.241	0.999
B36	B36	15/02/2022	y = 1.132x-3.625	0.996
B37	B37	04/02/2022	y = 1.157x+2.640	0.999
B38	B38	15/02/2022	y = 1.1432x-2.720	0.999
B39	B39	07/02/2022	y = 1.256x-7.614	1.000
B40	B40	15/02/2022	y = 1.175x-4.385	0.998
B41	B41	07/02/2022	y = 1.133x-1.951	0.998
B42	B42	04/02/2022	y = 1.127x-1.985	1.000
B43	B43	16/02/2022	y = 1.089x+0.223	0.996
B44	B44	03/02/2022	y = 1.339x-11.636	0.997
R01	R01	02/02/2022	y = 1.196x-5.960	0.996
R02	R02	09/02/2022	y = 1.175x-5.572	1.000
R03	R03	02/02/2022	y = 1.187x-6.283	0.995
R04	R04	07/02/2022	y = 1.100x-1.352	0.997
R05	R05	09/02/2022	y = 1.238x-8.500	0.997
R06	R06	01/02/2022	y = 1.328x-11.118	0.996
R07	R07	07/02/2022	y = 1.039x+1.507	0.995
R08	R08	04/02/2022	y = 1.141x-3.942	0.997
R09	R09	01/02/2022	y = 1.192x-5.710	0.997
R10	R10	09/02/2022	y = 1.194x-5.807	1.000
R11	R11	01/02/2022	y = 1.054x+0.098	0.996
R12	R12	04/02/2022	y = 1.171x-5.349	0.996
R13	R13	04/02/2022	y = 1.114x-1.755	0.999
R14	R14	07/02/2022	y = 1.100x-0.965	0.997
R15	R15	14/02/2022	y = 1.047x+1.073	0.995
R16	R16	09/02/2022	y = 1.129x-3.642	0.999
R17	R17	03/02/2022	y = 1.198x-5.739	1.000
R18	R18	02/02/2022	y = 1.268x-9.241	0.998
R19	R19	03/02/2022	y = 1.216x-5.626	0.999
R20	R20	01/02/2022	y = 1.197x-5.676	0.997

Calibrated by :





บริษัท เอส.พี.เอส. คอนซัลติ้ง เซอร์วิส จำกัด

S.P.S. CONSULTING SERVICE CO., LTD.

7 ซอยพหลโยธิน 24 ถนนพหลโยธิน แขวงจอมพล เขตจตุจักร กรุงเทพฯ 10900

7 Soi Phaholyothin 24, Phaholyothin Rd., Jompol, Chatuchak, Bangkok 10900

Tel : (662) 939-4370-72 Fax : (662) 513-4221 E-mail : sale@spscon.com, www.spscon.com

High Volume PM-10 Air Sampler Calibration Report

Calibration Method : Multipoint Orifice Flow Transfer Standard

Model : TE 5025A

S/N : 3095

Calibration Data

High Volume PM-10 Data		Calibration Data		
Recorder No.	Blower No.	Date	Actual Flowrate (ft ³ /min)	R ²
B01	B01	02/02/2022	y = 1.199x-0.729	0.999
B02	B02	04/02/2022	y = 1.047x+3.100	0.999
B03	B03	07/02/2022	y = 1.212x+3.555	0.997
B04	B04	09/02/2022	y = 1.314x-9.389	1.000
B05	B05	03/02/2022	y = 1.207x-5.472	0.995
B06	B06	04/02/2022	y = 1.260x-8.728	0.997
B07	B07	04/02/2022	y = 1.212x-5.353	0.996
B08	B08	09/02/2022	y = 1.285x-7.356	0.998
B09	B09	08/02/2022	y = 1.243x-6.277	1.000
B10	B10	07/02/2022	y = 1.285x-9.647	0.998
B11	B11	02/02/2022	y = 1.240x-6.135	0.995
B12	B12	01/02/2022	y = 1.285x-9.647	0.998
B13	B13	04/02/2022	y = 1.302x-9.419	0.996
B14	B14	07/02/2022	y = 1.199x+3.376	0.998
B15	B15	04/02/2022	y = 1.118x-0.993	0.995
B16	B16	04/02/2022	y = 1.190x-1.101	0.998
B17	B17	03/02/2022	y = 1.201x-2.953	0.998
B18	B18	07/02/2022	y = 1.143x-1.983	0.998
B19	B19	03/02/2022	y = 1.036x+1.865	0.998
B20	B20	03/02/2022	y = 1.201x-6.181	0.997
B21	B21	04/02/2022	y = 1.158x-0.828	0.998
B22	B22	04/02/2022	y = 1.290x-8.497	0.998
B23	B23	07/02/2022	y = 1.090x-0.542	1.000
B24	B24	01/02/2022	y = 1.218x-6.279	0.998
B25	B25	01/02/2022	y = 1.156x-3.313	0.997
B26	B26	07/02/2022	y = 1.135x+1.438	0.998
B27	B27	02/02/2022	y = 1.260x-8.474	0.998
B28	B28	04/02/2022	y = 1.090x-0.306	0.999
B29	B29	04/02/2022	y = 1.262x-8.639	1.000
B30	B30	03/02/2022	y = 1.219x-6.529	0.996
B31	B31	17/02/2022	y = 1.059x+0.716	0.997
B32	B32	16/02/2022	y = 1.154x-3.610	0.999
B33	B33	03/02/2022	y = 1.258x-8.776	0.999
B34	B34	16/02/2022	y = 1.123x+0.227	0.995

Calibrated



บริษัท เอส.พี.เอส. คอนซัลติ้ง เซอร์วิส จำกัด

S.P.S. CONSULTING SERVICE CO., LTD.

7 ซอยพหลโยธิน 24 ถนนพหลโยธิน แขวงจอมพล เขตจตุจักร กรุงเทพฯ 10900

7 Soi Phaholyothin 24, Phaholyothin Rd., Jompol, Chatuchak, Bangkok 10900

Tel : (662) 939-4370-72 Fax : (662) 513-4221 E-mail : sale@spscon.com, www.spscon.com

High Volume PM-10 Air Sampler Calibration Report

Calibration Method : Multipoint Orifice Flow Transfer Standard

Model : TE 5025A

S/N : 3095

Calibration Data

High Volume PM-10 Data		Calibration Data		
Recorder No.	Blower No.	Date	Actual Flowrate (ft ³ /min)	R ²
R01	R01	04/02/2022	y = 1.238x-7.598	0.995
R02	R02	11/02/2022	y = 1.161x-3.899	0.996
R03	R03	04/02/2022	y = 1.154x+2.827	0.998
R04	R04	06/02/2022	y = 1.116x-1.752	0.995
R05	R05	07/02/2022	y = 1.125x-2.487	0.995
R06	R06	10/02/2022	y = 1.321x-9.065	0.998
R07	R07	04/02/2022	y = 1.138x-1.986	0.996
R08	R08	03/02/2022	y = 1.160x-3.759	0.996
R09	R09	10/02/2022	y = 1.209x-6.918	0.995
R10	R10	04/02/2022	y = 1.114x-1.889	0.995
R11	R11	03/02/2022	y = 1.272x-7.597	1.000
R12	R12	03/02/2022	y = 1.153x-3.385	0.995
R13	R13	02/02/2022	y = 1.207x-4.913	0.996
R14	R14	01/02/2022	y = 1.183x-3.660	0.996
R15	R15	02/02/2022	y = 1.247x-7.741	0.999
R16	R16	02/02/2022	y = 1.238x-6.677	0.996
R17	R17	01/02/2022	y = 1.203x-5.310	0.998
R18	R18	04/02/2022	y = 1.148x-3.211	0.998
R19	R19	04/02/2022	y = 1.220x-6.839	0.997
R20	R20	03/02/2022	y = 1.161x-5.047	0.997

Calibrated by :

(Mr. Phakhinai Khongkomnerd)

Approved by :

(Mr. Peera Detudom)

**QUALITY CALIBRATION CO.,LTD.**

235 Petchkasem 63/2 Road, Laksong, Bangkae, Bangkok 10160

Tel (662) 421-5402, (662) 444-0152-3, Fax (662) 809-4584

www.qcalibration.comNSC-TISI-TISI7025
CALIBRATION 0049

CERTIFICATE No : 21M3169

REFERENCE No : 60627-5

PAGE : 1 OF 2

Certificate of Calibration

EQUIPMENT : DIGITAL BALANCE

MANUFACTURER : METTLER TOLEDO

MODEL : XS105DU

SERIAL No : 1126422905

ID No : BA 05/50

CONDITION AS RECEIVED : USED ITEM

SUBMITTED BY : S.P.S. CONSULTING SERVICE CO., LTD.
7 SOI PHAHOLYOTHIN 24, PHAHOLYOTHIN RD.,
JOMPOL, CHATUCHAK, BANGKOK 10900

CALIBRATED BY : ATSAWIN Y.

CALIBRATION DATE : 19-Mar-21

APPROVED BY : 

ISSUED DATE : 20-Mar-21

RECEIVED DATE : 19-Mar-21

THIS CERTIFICATE MAY NOT BE REPRODUCED OTHER THAN IN FULL EXCEPT WITH THE PRIOR WRITTEN APPROVAL OF
QUALITY CALIBRATION CO., LTD.

**QUALITY CALIBRATION CO.,LTD.**

235 Petchkasem 63/2 Road, Laksong, Bangkai, Bangkok 10160

Tel (662) 421-5402, (662) 444-0152-3, Fax (662) 809-4584

www.qcalibration.com

CERTIFICATE No : 21M3169

PAGE : 2 OF 2

Calibration Report

EQUIPMENT : DIGITAL BALANCE MODEL : XS105DU
MANUFACTURER : METTLER TOLEDO S/N : 1126422905
ID No : BA 05/50 RECEIVED DATE : 19-Mar-21
AIR PRESSURE : 1009mbar \pm 1mbar CALIBRATION DATE : 19-Mar-21
AMBIENT TEMPERATURE : 24°C \pm 1°C RELATIVE HUMIDITY : 52 %RH \pm 10 % RH

CONDITION OF THIS RESULTS OF CALIBRATION

1. THIS INSTRUMENT WAS CALIBRATED BY ACCORDING TO UKAS LAB 14 EDITION 6:2019 BY USING KNOWN WEIGHT STANDARD WEIGHT. THE BALANCE WAS ADJUSTED USING WEIGHT OF QUALITY CALIBRATION TO ADJUST. THE BALANCE HAS NO ZERO TRACKING FUNCTION. REPEATABILITY WAS MEASURED BY USING 10 REPEATED MEASUREMENTS. LINEARITY WAS MEASURED COVERING 10 POINTS, EVENLY SPREAD OVER THE RANGE. THE INSTRUMENT WAS SET ZERO BEFORE PERFORMING THE LINEARITY TEST. OFF-CENTER LOADING WAS MEASURED BY USING STANDARD WEIGHTS PLACED ON THE PAN AND MOVED TO VARIOUS POSITIONS ON THE PAN.

2. REFERENCE STANDARD INSTRUMENTS :-

INSTRUMENT	MODEL	SERIAL No	CERTIFICATE No	DUE DATE
1) STANDARD WEIGHT SET	E2	QK-I-151	C02210415	09-Feb-23
2) STANDARD WEIGHT	E2	15843	C02210419	10-Feb-23
3) STANDARD WEIGHT	E2	QK-I-349	M2103235S	26-Mar-23

3. THIS RESULT WAS FOUND ACCURATE AS SHOWN ON DATE AND PLACE OF CALIBRATION ONLY.

4. THIS RESULT EXCLUDE LONG TERM STABILITY OF THE UNIT UNDER CALIBRATION.

5. THIS CERTIFICATE IS TRACEABLE TO THE INTERNATIONAL SYSTEM OF UNIT MAINTAINED AT:-

- NATIONAL INSTITUTE OF METROLOGY (THAILAND) THROUGH CENTRAL BUREAU OF WEIGHTS&MEASURES

RESULT OF CALIBRATION :- WITHOUT ADJUSTMENT


1. ZERO SETTING FUNCTION : NORMAL

2. TARE FUNCTION : NORMAL

3. REPEATABILITY OF READING AT 100 g WAS 0.000055 g

4. DEPARTURE FROM NOMINAL VALUE/ LINEARITY

NOMINAL VALUE (g)	BALANCE READING (g)	CORRECTION (g)	UNCERTAINTY (\pm g)
0.00	0.00000	0.00000	0.000066
0.02	0.01998	0.00002	0.000066
0.10	0.10001	-0.00001	0.000066
0.20	0.20001	-0.00001	0.000067
0.50	0.49996	0.00004	0.000065
1.00	0.99997	0.00003	0.000066
2.00	2.00000	0.00000	0.000067
5.00	5.00002	-0.00002	0.000068
10.00	10.00003	-0.00003	0.000070
20.00	20.00000	0.00000	0.000075
50.00	50.00000	0.00000	0.00013
100.00	100.0001	-0.0001	0.00019
120.00	120.0001	-0.0001	0.00022

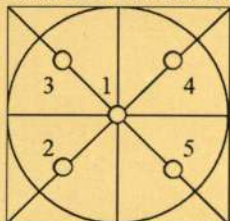5. OFF CENTER LOADING ERROR

POINT	READING (g)
1	50.0000
2	50.0000
3	50.0000
4	50.0000
5	50.0000
OFF-CENTER LOADING	0.0000

NOTE: THIS CALIBRATION WAS CARRIED OUT AT THE CUSTOMER'S PLACE AT LABORATORY AREA

THE REPORTED UNCERTAINTY OF MEASUREMENT WAS BASED ON A STANDARD UNCERTAINTY MULTIPLIED BY A COVERAGE FACTOR $k=2$, PROVIDING A LEVEL OF CONFIDENCE APPROXIMATELY 95%.

END OF CALIBRATION REPORT

THAILAND INSTITUTE OF SCIENTIFIC AND TECHNOLOGICAL RESEARCH (TISTR)

Request No. 21-64/0528

MTC No. EEL. BP. 17/0564

CALIBRATION CERTIFICATE

Submitted by : S.P.S. Consulting Services Service Co.,Ltd.

Address : 7 Soi Phaholyothin 24, Phaholyothin Road, Jompol, Chatuchak, Bangkok 10900.

Calibrated at : Electrical and Electronic Standards Laboratory, Industrial Metrology and Testing Service Centre.
: Soi 1C, Bangpoo Industrial Estate, Sukhumvit Rd., Muang, Samutprakan 10280.

Instrument Calibrated :

Description : Sound Calibrator

Manufacturer : ACO

Model : 2127

Serial No. : 130006

Ambient Environment

Temperature : $(23 \pm 3) ^\circ\text{C}$

Relative Humidity : $(50 \pm 15) \%$

Ambient Pressure : $(101.325 \pm 1.500) \text{ kPa}$

Standards used : 1. Digital Function Synthesizer NF Electronic DF-193A S/N 122037.
2. Measuring Amplifier Bruel&Kjaer 2636 S/N 1537484.
3. Programmable Attenuator Tamagawa TPA-303A S/N OF 2214.
4. Digital Multimeter Agilent 34401A S/N MY44005560.
5. Pressure Transmitter Vaisala PTB202AD S/N T0650001.
6. Audio Analyzer Keithley 2015-P S/N 4106495.
7. Condenser Microphone Bruel&Kjaer 4180 S/N 2889871.

Calibration Procedure: CP-102-04 based on IEC 60942-2003; The sound pressure level generated by sound calibrator under test shall be measured by standard microphone using an insert voltage technique.

This instrument has been calibrated against standards maintained at Electrical and Electronic Standards Laboratory (EEL), which are traceable to the International System of Units through the National Institute of Metrology (Thailand).

The information on actual reading is attached herewith and the uncertainty limits quoted refer to the measured values only.

Date of Receipt : 6 May 2021

Date of Calibration : 15 May 2021

1 / 2 ✓

The results relate only to the items tested/calibrated or value assigned.

Advertising the Report/Certificate and publicity of the results except in full are prohibited unless written permission is obtained from the governor of TISTR.

FM.BL.MTC.002 Rev.4

Head Office

35 Mu 3 Tambon Khlone Ha, Amphoe Khlong Luang,
Changwat Pathumthani 12120, Thailand

Tel. (66) 0 2577 9000

Fax. (66) 0 2577 9009

E-mail : rumpai@tistr.or.th Website:www.tistr.or.th

Office/Laboratory

Soi 1C, Bangpoo Industrial Estate, Sukhumvit Road,
Amphoe Muang, Changwat Samutprakan 10280, Thailand

Tel. (66) 0 2323 1672-80 ext. 115, 116

Fax. (66) 0 2323 9165

E-mail : mtc@tistr.or.th

Office

196 Phahonyothin Road, Chatuchak, Bangkok 10900,
Thailand

Tel. (66) 0 2579 1121-30 ext. 5219, 5225, 5217

Fax. (66) 0 2579 8592

E-mail : sumalee@tistr.or.th

THAILAND INSTITUTE OF SCIENTIFIC AND TECHNOLOGICAL RESEARCH (TISTR)

Request No. 21-64/0528

MTC No. EEL. BP. 17/0564

The reported expanded uncertainty is based upon a standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%.

Nominal Output of Unit Under Test = 94 dB re 20 μ Pa at 1000 Hz

Acoustic Output in dB re 20 μ Pa, Corrected to Reference Conditions: 101.325 kPa, 23.0 °C and 50 %RH.

1. Sound Pressure Level

Standard Microphone Type	Measured Sound Pressure Level (dB)	Deviated value (dB)	Uncertainty (dB)	Tolerance limit IEC60942:2003 Class 1
1/2 inch Brüel&Kjaer 4180	93.96	-0.04	± 0.10	± 0.40 dB

2. Frequency

Standard Microphone Type	Measured Frequency (Hz)	Deviated value (Hz)	Uncertainty (Hz)	Tolerance limit IEC60942:2003 Class 1
1/2 inch Brüel&Kjaer 4180	999.9	-0.1	± 1.5	$\pm 1.0\%$

3. Total Distortion

Standard Microphone Type	Measured Total Distortion (%)	Uncertainty (%)	Tolerance limit IEC60942:2003 Class 1
1/2 inch Brüel&Kjaer 4180	1.26	± 0.50	$\pm 3.0\%$

Note : 1. No adjustment.

2. The calibrator pressure correction was not included.

3. The microphone volume correction was not included.

Calibrated by :



Acting Director

Electrical and Electronic Standards Laboratory
Industrial Metrology and Testing Service Centre

Date of Calibration : 15 May 2021

Date of Issue : 18 May 2021

Ref : 2011264050601894002

End of Certificate

2 / 2

The results relate only to the items tested/calibrated or value assigned.

Advertising the Report/Certificate and publicity of the results except in full are prohibited unless written permission is obtained from the governor of TISTR.

FM.BLMTC.002 Rev.4

Head Office

35 Mu 3 Tambon Khlong Ha, Amphoe Khlong Luang,
Changwat Pathumthani 12120, Thailand
Tel. (66) 0 2577 9000
Fax. (66) 0 2577 9009
E-mail : rumpai@tistr.or.th Website:www.tistr.or.th

Office/Laboratory

Soi 1C, Bangpoo Industrial Estate, Sukhumvit Road,
Amphoe Muang, Changwat Samutprakan 10280, Thailand
Tel. (66) 0 2323 1672-80 ext. 115, 116
Fax. (66) 0 2323 9165
E-mail : mtc@tistr.or.th

Office

196 Phahonyothin Road, Chatuchak, Bangkok 10900,
Thailand
Tel. (66) 0 2579 1121-30 ext. 5219, 5225, 5217
Fax. (66) 0 2579 8592
E-mail : sumalee@tistr.or.th



บริษัท เอส.พี.เอส. คอนซัลติ้ง เซอร์วิส จำกัด
S.P.S. CONSULTING SERVICE CO., LTD.
7 ซอยพหลโยธิน 24 ถนนพหลโยธิน แขวงจอมพล เขตจตุจักร กรุงเทพฯ 10900
7 Soi Phaholyothin 24, Phaholyothin Rd., Jompol. Chatuchak, Bangkok 10900
Tel : (662) 939-4370-72, Fax : (662) 513-4221 E-mail : sale@spscon.com., www.spscon.com

Noise B_010/22

Sound Level Meter Calibration Report


Acoustic Calibrator Data

Brand	ACO	Number	AC 03/56
Model	2127	Serial No.	130006
Calibration Range	94 dB, 1000 Hz	Last Calibration	15 May 2021
		Due Date	15 May 2022

Calibration Data

Sound Level Meter Data				Calibration Data		
SLM No.	Brand	Model	Serial No.	Date	Actual Reading [dB]	
					Before Adjustment	After Adjustment
ACO-B06	ACO	6236	00142003	14 March 2022	94.0	94.0
ACO-B07	ACO	6236	00142004	14 March 2022	94.0	94.0
ACO-B17	ACO	6236	00172042	14 March 2022	93.9	94.0
ACO-B21	ACO	6236	00172059	14 March 2022	93.9	94.0
ACO-B40	ACO	6236	00192031	14 March 2022	94.0	94.0
Acoustic Certified Value : Thailand Institute of Scientific and Technological Research (TISTR)					93.96 ± 0.40 dB	

Calibrated by :

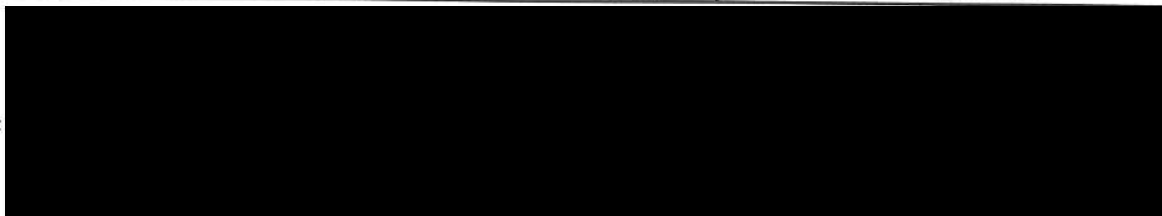



CERT NO.: C22/0066B

PAGE : 1 OF 3

**ISOCAL TECHNOLOGY CO.,LTD.
INDUSTRIAL INSTRUMENT CALIBRATION CENTER**

170/405 Moo 3 Serithai Rd., Kannayao Kannayao Bangkok 10230

Tel. 0-2906-3040-1 Fax. 0-2919-9948

Certificate of Calibration

EQUIPMENT :	PH METER
MODEL :	HI98190
SERIAL NO. :	04260035101
ID NO. :	B01
MANUFACTURER :	HANNA
MADE IN :	ROMANIA
SUBMITTED BY :	S.P.S. CONSULTING SERVICE CO.,LTD. 7 SOI PHAHOLYOTHIN 24, PHAHOLYOTHIN RD.,JOMPOL, CHATUCHAK , BANGKOK
AMBIENT TEMPERATURE :	(23 ± 2) °C
RELATIVE HUMIDITY :	(50 ± 15) %
CALIBRATED BY:	
APPROVED BY :	
ISSUE DATE :	23-Mar-2022

THIS CERTIFICATE MAY NOT BE REPRODUCED OTHER THAN IN FULL, EXCEPT WITH THE PRIOR
WRITTEN APPROVAL OF THE HEAD OF THE INDUSTRIAL INSTRUMENTS CALIBRATION CENTER.



ISOCAL TECHNOLOGY CO., LTD.

CALIBRATION REPORT

CERT NO.: C22/0066B

PAGE : 2 OF 3

EQUIPMENT : PH METER
MANUFACTURER : HANNA
MODEL : HI98190
SERIAL NO. : 04260035101
ID NO. : B01
CALIBRATION DATE : 19-Mar-2022
RECEIVED DATE : 17-Mar-2022
PROCEDURE USED :

CALIBRATION WERE CONDUCTED USING IN-HOUSE CALIBRATION PROCEDURE WI-18-22 ACCORDING TO COMPARISON WITH PH SOLUTION STANDARD.

CONDITION OF THIS RESULT OF CALIBRATION

1. THIS RESULT OF CALIBRATION WAS FOUND ACCURATE AS SHOWN ON DATE AND PLACE OF CALIBRATION ONLY.
2. THIS CERTIFICATE IS TRACEBLE TO THE INTERNATIONAL OF THIS RESULT OF CALIBRATION.
3. REFERENCE STANDARDS INSTRUMENTS :-

PH SOLUTION MODEL PH 4.01 SERIAL NO. 1.09435.1000 CERT. NO HC02910835 DATE 24-APR-2022

- NATIONAL INSTITUTE OF STANDARDS AND TECHNOLOGY(NIST), U.S.A
-PHYSIKALISCH-TECHNISCHE BUNDESANSTALT (PTB),GERMANY.
THROUGH SUPELCO CO., LTD.

PH SOLUTION MODEL PH 7.01 SERIAL NO. 1.09439.1000 CERT. NO HC02387439 DUE DATE 18-APR-2022

- NATIONAL INSTITUTE OF STANDARDS AND TECHNOLOGY(NIST), U.S.A
-PHYSIKALISCH-TECHNISCHE BUNDESANSTALT (PTB),GERMANY.
THROUGH SUPELCO CO., LTD.

PH SOLUTION MODEL PH 10.01 SERIAL NO. 1.09438.1000 CERT. NO HC01501438 DUE DATE 27-MAR-2022

- NATIONAL INSTITUTE OF STANDARDS AND TECHNOLOGY(NIST), U.S.A
-PHYSIKALISCH-TECHNISCHE BUNDESANSTALT (PTB),GERMANY.
THROUGH SUPELCO CO., LTD.



ISOCAL TECHNOLOGY CO., LTD.

CALIBRATION REPORT

CERT.NO.: C22/0066B

PAGE : 3 OF 3

RESULT OF CALIBRATION: ADJUSTMENT (YES)

FUNCTION: PH MEASUREMENT WITH SOLUTION @ 25 °C

SCALE RANGE : 4.01 pH TO 10.01 pH

RESOLUTION: 0.01 pH

STANDARD VALUE (pH)	UUC READING (pH)	ERROR (pH)	UNCERTAINTY (pH)
4.01	3.97	-0.04	0.012
7.01	7.03	0.02	0.012
10.01	10.04	0.03	0.012

THE REPORTED UNCERTAINTY OF MEASUREMENT WAS BASED ON STANDARD UNCERTAINTY MULTIPLIED BY A COVERAGE FACTOR $k=2$, PROVIDING A LEVEL OF CONFIDENCE OF APPROXIMATELY 95%

UUC = UNIT UNDER CALIBRATE

- oOo -



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CALIBRATION AND TESTING EQUIPMENT SERVICES

534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250

TEL. 0-2717-3000-24 FAX. 0-2719-9484

Cert.No.: 21CH52

Page.: 1 of 2

Certificate of Calibration

Equipment : Turbidity Meter
Manufacturer : Thermo Scientific
Model : EUTECH TN-100
Serial No. : 2882693
ID. No. : TB 04/62
Condition As-Received: Used Item
Received Date : 15 January 2021
Calibration Date : 18 January 2021
Reference : 2101-0429WN-2
Submitted by : S.P.S. Consulting Service Co.,Ltd.
7 Soi Phaholyothin 24, Phaholyothin Rd.,
Jompol, Chatuchak, Bangkok 10900
Ambient Temperature : (25 \pm 2.5) °C
Relative Humidity : (50 \pm 20) %
Calibration Procedure : In - house method : CP-CH11
based on direct measurement by
using Formazin standard solution

Calibrated by :

Approved by :

(/) Malee Butkruea
() Saithip Meangmai
() Warakorn Lerngagtrakul

Issue Date : 26 January 2021

The Uncertainties are for a confidence probability of approximately 95%.

This certificate may not be reproduced other than in full, except with the prior written
approval of the head of Calibration and Testing Equipment Services.

A 0005766



Cert.No. : 21CH52

Page. : 2 of 2

Condition of this calibration result

1. Reference Standard Instruments :

This certification is traceable to the International System of unit through the reference standards laboratory of Industrial Calibration Center, Technology Promotion Association (Thailand-Japan).

<u>Instruments</u>	<u>Serial No.</u>	<u>ID No.</u>	<u>Certificate No.</u>	<u>Due date</u>
1) Thermo-Hygrograph	1103328	130EC010	20H1607	02 July 2021
2) Electronic Balance	1126143764	140RC004	20MM595	27 Sep 2021

2. Standard Material : The Formazin suspension has been prepared gravimetric from

<u>Material</u>	<u>Manufacturer</u>	<u>Lot No.</u>	<u>Assay</u>
1) Hexamethylenetetramine	HIMEDIA	0000343342	99.5%
2) Hydrazinium Sulfate	HIMEDIA	0000332928	99.2%

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration result

**Performing three - Formazin suspension standard curve by using 20,100,800 NTU
Turbidity Meter Serial Number : 2882693**

Standard Formazine suspension (NTU)	UUC* Reading (NTU)	Uncertainty of Measurement (± NTU)	Coverage Factor <i>k</i>
20	20.0	0.38	2.00
40	40.2	0.40	2.00
100	99.6	0.72	2.00
400	457	1.5	2.07

Remark

- UUC* = Unit Under Calibration
- NTU = Nephelometric Turbidity Units

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

malu

a 1037061



CERTIFICATE No : 21M3167

REFERENCE No : 60627-3

PAGE : 1 OF 2

Certificate of Calibration

EQUIPMENT : DIGITAL BALANCE

MANUFACTURER : SARTORIUS

MODEL : BSA224S-CW

SERIAL No : 36591843

ID No : BA 09/61

CONDITION AS RECEIVED : USED ITEM

SUBMITTED BY : S.P.S. CONSULTING SERVICE CO., LTD.
7 SOI PHAHOLYOTHIN 24, PHAHOLYOTHIN RD.,
JOMPOL, CHATUCHAK, BANGKOK 10900

CALIBRATED BY : ATSAWIN Y.

CALIBRATION DATE : 19-Mar-21

APPROVED BY : 

ISSUED DATE : 20-Mar-21

RECEIVED DATE : 19-Mar-21



CERTIFICATE No : 21M3167

PAGE : 2 OF 2

Calibration Report

EQUIPMENT : DIGITAL BALANCE MODEL : BSA224S-CW
MANUFACTURER : SARTORIUS S/N : 36591843
ID No : BA 09/61 RECEIVED DATE : 19-Mar-21
AIR PRESSURE : 1009mbar \pm 1mbar CALIBRATION DATE : 19-Mar-21
AMBIENT TEMPERATURE : 24°C \pm 1°C RELATIVE HUMIDITY : 52 %RH \pm 10 % RH

CONDITION OF THIS RESULTS OF CALIBRATION

1. THIS INSTRUMENT WAS CALIBRATED BY ACCORDING TO UKAS LAB 14 EDITION 6:2019 BY USING KNOWN WEIGHT STANDARD WEIGHT. THE BALANCE WAS ADJUSTED USING WEIGHT OF QUALITY CALIBRATION TO ADJUST. THE BALANCE HAS NO ZERO TRACKING FUNCTION. REPEATABILITY WAS MEASURED BY USING 10 REPEATED MEASUREMENTS. LINEARITY WAS MEASURED COVERING 10 POINTS, EVENLY SPREAD OVER THE RANGE. THE INSTRUMENT WAS SET ZERO BEFORE PERFORMING THE LINEARITY TEST. OFF-CENTER LOADING WAS MEASURED BY USING STANDARD WEIGHTS PLACED ON THE PAN AND MOVED TO VARIOUS POSITIONS ON THE PAN.

2. REFERENCE STANDARD INSTRUMENTS :-

INSTRUMENT	MODEL	SERIAL No	CERTIFICATE No	DUE DATE
1) STANDARD WEIGHT SET	E2	QK-I-151	C02210415	09-Feb-23
2) STANDARD WEIGHT	E2	15843	C02210419	10-Feb-23
3) STANDARD WEIGHT	E2	QK-I-349	M2103235S	26-Mar-23

3. THIS RESULT WAS FOUND ACCURATE AS SHOWN ON DATE AND PLACE OF CALIBRATION ONLY.

4. THIS RESULT EXCLUDE LONG TERM STABILITY OF THE UNIT UNDER CALIBRATION.

5. THIS CERTIFICATE IS TRACEABLE TO THE INTERNATIONAL SYSTEM OF UNIT MAINTAINED AT:-

- NATIONAL INSTITUTE OF METROLOGY (THAILAND) THROUGH CENTRAL BUREAU OF WEIGHTS&MEASURES

RESULT OF CALIBRATION :- WITHOUT ADJUSTMENT

1. ZERO SETTING FUNCTION : NORMAL


2. TARE FUNCTION : NORMAL

3. REPEATABILITY OF READING AT 200 g WAS 0.000045 g

4. DEPARTURE FROM NOMINAL VALUE/ LINEARITY

NOMINAL VALUE (g)	BALANCE READING (g)	CORRECTION (g)	UNCERTAINTY (\pm g)
0.0	0.0000	0.0000	0.000075
0.1	0.1000	0.0000	0.000075
0.2	0.2000	0.0000	0.000076
0.5	0.5000	0.0000	0.000076
1.0	1.0000	0.0000	0.000077
2.0	2.0000	0.0000	0.000077
5.0	5.0000	0.0000	0.000079
10.0	10.0000	0.0000	0.000082
20.0	20.0000	0.0000	0.000086
50.0	50.0000	0.0000	0.00013
100.0	100.0001	-0.0001	0.00019
200.0	199.9997	0.0003	0.00032

5. OFF CENTER LOADING ERROR

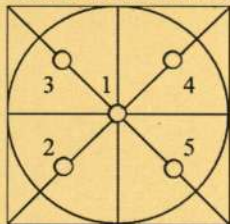

POINT	READING (g)
1	100.0000
2	100.0000
3	100.0001
4	100.0000
5	99.9999
OFF-CENTER LOADING	0.0001


NOTE: THIS CALIBRATION WAS CARRIED OUT AT THE CUSTOMER'S PLACE AT LABORATORY AREA

THE REPORTED UNCERTAINTY OF MEASUREMENT WAS BASED ON A STANDARD UNCERTAINTY MULTIPLIED BY A COVERAGE FACTOR $k=2$, PROVIDING A LEVEL OF CONFIDENCE APPROXIMATELY 95%.

END OF CALIBRATION REPORT

Lambda UV Preventive Maintenance (PM)

Company Name:	S.P.S. CONSULTING SERVICE CO., LTD.		
Address:	7, Soi Phaholyothin24, Ladyao, Jatujak, Bangkok		
User Name:	K. Benjawan	WO Number:	WO-01550999
Telephone Number:	086-141-2523	PM Number:	6 of 6 P
Customer Support Engineer:	K. Anon	Certificate Number:	UV2004-2022
Date PM Performed: (DD-MMM-YYYY)	25-Jan-2022	Next PM Due Date: (DD-MMM-YYYY)	25-Jul-2022

Part Number	Release	Publication Date	
09370504	B	March 2013	

Scope

The purpose of this PM is to ensure the continued functionality of the PerkinElmer Lambda UV/Vis Spectrophotometer by inspecting and replacing any worn or damaged parts. This service should only be performed by a trained representative of PerkinElmer. The customer should save their method before the PM begins.

General Instructions:

The customer must provide the engineer operational data to demonstrate recent instrument performance prior to starting the PM. Always check with the customer before making any changes that may affect the customer's analysis should be signed by an authorized PerkinElmer and customer representative and left with the customer. Update the PM sticker and instrument logbook as required.

Copyright Information

This document contains proprietary information that is protected by copyright. All rights are reserved. No part of this publication may be reproduced in any form whatsoever or translated into any language without the prior, written permission of PerkinElmer, Inc. Copyright © 2009 PerkinElmer, Inc.

Trademarks

Registered names, trademarks, etc. used in this document, even when not specifically marked as such, are protected by law. PerkinElmer is a registered trademark of PerkinElmer, Inc. All other trademarks and registered trademarks not owned by PerkinElmer, Inc. or its subsidiaries that are depicted herein are the property of their respective owners. Except as specifically set forth in its terms and conditions of sale, PerkinElmer makes no Warranty of any kind with regard to this document, including, but not limited to, the implied warranties of merchantability and fitness for a particular purpose. PerkinElmer shall not be liable for incidental or consequential damages in connection with the furnishing or use of this document.

Component List

Component Specific Model	Serial #	Software Version		Configuration Notes
Lambda 25	501S14123010	6.2.0.0741	STD	1.27
NA	NA	NA	NA	NA

Parts Lists

Parts Included with the PM				
Part Number (if applicable)	Description	Quantity	Serial Number	Expiration Date (MM/YY)
B250 0099	Stray Light standard			
	Nal cell	1	1943	Jan-22
	NaNO2 cell	1	2963	
	KCl cell	1	31030	
	H2O	1	71497	
B050 7805	Secondary Standards for calibration of wavelength and photometric accuracy or use NBS/NIST 390 standards			
	Gray Glass G1	1	2926	Jan-22
	Gray Glass G2	1	3501	
	Gray Glass G3	1	2552	
	Holmium Glass	1	1085	

Additional Tools Required for PM					
Part Number (if applicable)	Description	Quantity	Serial #		Remark
-	-	-	-	-	-
-	-	-	-	-	-
-	-	-	-	-	-
Additional Reagents and Standards Required for PM					
Part Number (if applicable)	Description	Quantity	Batch/Lot #		Expiration Date (MM/YY)
-	-	-	-	-	-
-	-	-	-	-	-
-	-	-	-	-	-

Procedure Checklist

Use (✓) to check off those steps in the checklist that have been completed.

1. General:

- ☒ Review the instrument performance with the customer and document any recent problems.
- ☒ Inspect the customer log book and make any appropriate PM entries.
- ☒ Perform general inspection of system for cleanliness.

2. Optical checks:

- ☒ Lamp Alignment/Energy
- ☒ Sample Compartment Windows/Monochromator
- ☒ Mirror and Grating Alignment
- ☒ Cell Holder Alignment

3. Mechanical:

- ☒ Physical inspection – Please write any comments in the additional comments section.
- ☒ Grating Drive Mechanism.
- ☒ Lamp Change Mechanism.
- ☐ Slit Drive Manual Servo.

4. Test:

Refer to Appendix A for the specifications of the instrument being tested.

- ☒ D2 Wavelength accuracy

	Actual Value	Specification
Accuracy at 656.1 nm	656.16	± 0.1

☒ Holmium Oxide wavelength accuracy

Filter ID #		1085		
Test	Calibration Value	Actual Value	Deviation	Specification
279.3 nm	279.3	279.39	-0.09	± 0.5
360.8 nm	360.9	360.93	-0.03	± 0.5
459.9 nm	460.0	460.07	-0.07	± 0.5
536.4 nm	536.2	536.40	-0.20	± 0.5

☒ Scattered Light.

Test	Filter ID #	Result	Specification
NaI @ 220 nm	1943	0.0133	< 0.02 %T
NaNO ₂ @ 340 nm	2963	-0.1296	< 0.02 %T
NaNO ₂ @ 370 nm	2963	-0.0002	< 0.02 %T
KCl @ 200 nm	31030	2.4808	≥ 2 A

☒ Baseline Flatness.

Corrected Baseline	Specification
0.000163	± 0.001 A

☒ Noise Test @ 500 nm.

Actual Value	Specification
0.0000240	± 0.00008 A

☒ Photometric Accuracy.

Filter 1 ID #		2926		
Test	Calibrated Value	Actual Value	Deviation	Specification
440 nm	0.3483	0.3493	-0.0010	± 0.006 A
546 nm	0.3029	0.3046	-0.0017	± 0.006 A
635 nm	0.3200	0.3232	-0.0032	± 0.006 A
Filter 2 ID #		3501		
Test	Calibrated Value	Actual Value	Deviation	Specification
440 nm	1.001	1.0024	-0.0014	± 0.006 A
546 nm	0.9797	0.9813	-0.0016	± 0.006 A
635 nm	1.0285	1.0325	-0.0040	± 0.006 A
Filter 3 ID #		2552		
Test	Calibrated Value	Actual Value	Deviation	Specification
440 nm	0.489	0.4935	-0.0045	± 0.006 A
546 nm	0.4582	0.4595	-0.0013	± 0.006 A
635 nm	0.5046	0.5075	-0.0029	± 0.006 A

5. Accessory (where applicable):

- ☐ Integrating Sphere
- ☐ Reflecting Attachment
- ☐ Cell Changer
- ☐ Sipper
- ☐ Auto Sampler


6. Review:

- ☒ Review with the customer PM work performed.
- ☒ Review with the customer routine maintenance procedures.
- ☒ Discuss recommended customer-supplied materials to have on hand
- ☒ Attach PM sticker.
- ☒ Update Logbook.

Additional Comments

Additional Comments Regarding the PM

Review

<p><i>The preventive maintenance checks and if applicable performance tests for Lambda UV have been completed.</i></p>	
<p>This Lambda UV Passes <input checked="" type="checkbox"/> Fails <input type="checkbox"/> the preventive maintenance.</p>	
<p>Review of Preventive Maintenance:</p>	
<p>Authorized PerkinElmer Representative:</p> <p>Anon Leenthawonkit </p>	<p>Date:</p> <p>25-Jan-2022 (DD-MM-YYYY)</p>
<p>Authorized Customer Representative:</p>	<p>Date:</p> <p>25-Jan-2022 (DD-MM-YYYY)</p>



MAINTENANCE AND TEST CERTIFICATE MODEL

OPTIMA 5300DV

Customer : <u>S.P.S.Consulting Service Co.,Ltd</u>	Date Tested: <u>January 12, 2022</u>	
	Recommendation Recertification	
Address : <u>7 Soi Phaholyothin 24</u>	Period <u>6</u> Months	
<u>Paholyothin Road</u>	Recertification Due: <u>July 12, 2022</u>	
<u>Jompol Chatuchak, Bangkok 1090</u>	Date Last Certified: <u>July 14, 2021</u>	
User Name: <u>K.Phenpha Viphasathawat</u>	Visit Number: <u>2 of 2</u>	
Phone: <u>083-9269252</u>	PerkinElmer Phone: <u>02-719-6420 ext 206</u>	
Fax: <u>02-513-4221</u>	PerkinElmer Fax: <u>02-318-5597</u>	

CONFIGURATION TESTED		ACCESSORIES/COMPONENT NOT INCLUDED
MODEL	SERIAL NUMBER	
<u>OPTIMA 5300DV</u>	<u>077C7042401</u>	
TESTED EQUIPMENT	CALIBRATION NUMBER	EXPIRATION
<u>IPV Methods</u>		
TEST STANDARD USED	PART NUMBER	EXPIRATION DATE
<u>Multielement Standard</u>	<u>N069-1579</u>	<u>August 30, 2022</u>
<u>Wavecal Solution</u>	<u>N058-2152</u>	<u>January 30, 2022</u>
<u>VIS Wavecal solution</u>	<u>N930-2946</u>	<u>June 30, 2022</u>
<u>Instrument Cal. STD4</u>	<u>N930-0221</u>	<u>August 30, 2022</u>
CUSTOMER SUPPLIED	COMMENTS	CUSTOMER INITIALS
<u>2 % HNO3</u>		
<u>10 % HNO3</u>		



MAINTENANCE AND TEST CERTIFICATE MODEL

OPTIMA 5300DV

SERIAL NUMBER 077C8011701

DATE TESTED January 12, 2022

1. MECHANICAL CHECKS

- | | |
|--|-----------------------------|
| A. Inspect and clean all fans and filters. | <input type="checkbox"/> OK |
| B. Inspect and replace as necessary, all torch components including the RF coil. | <input type="checkbox"/> OK |
| C. Inspect all tubing for sign of clacking or leaking. | <input type="checkbox"/> OK |
| D. Adjust water and gas pressure regulator settings. | <input type="checkbox"/> OK |
| E. Inspect and leak check pneumatics drawers. | <input type="checkbox"/> OK |
| F. Clean the exterior of the instrument. | <input type="checkbox"/> OK |

2. OPTICAL CHECKS

- | | |
|---|-----------------------------|
| A. Inspect and clean all optical components. | <input type="checkbox"/> OK |
| B. As required, check and replace all purgefilters. | <input type="checkbox"/> OK |
| C. Recheck optical alignment. | <input type="checkbox"/> OK |

3. COOLING SYSTEM CHECKS

- | | |
|---|------------------------------|
| A. Perform preventive maintenance on chiller. | <input type="checkbox"/> OK |
| B. Flush out the chiller every year. | <input type="checkbox"/> N/A |

4. PERFORMANCE CHECKS

- | | |
|----------------------------|-----------------------------|
| A. Torch View Alignment. | <input type="checkbox"/> OK |
| B. Wavelength Calibration. | <input type="checkbox"/> OK |



MAINTENANCE AND TEST CERTIFICATE MODEL

OPTIMA 5300DV

SERIAL NUMBER : 077C8011701

DATE TESTED : January 12, 2022

PARAMETER	SPECIFICATION			FINAL VALUE	
Spectral Resolution : UV	As	193.696 nm	≤ 0.007	0.00554	
	Ni	231.604 nm	≤ 0.008	0.00725	
	Ni	341.476 nm	≤ 0.012	0.00752	
Spectral Resolution : VIS	La	408.672 nm	≤ 0.020	0.01616	
	Ba	455.403 nm	≤ 0.025	0.02416	
Precision					
	As	193.656 nm	% RSD < 1.0	0.34	%
	Zn	213.856 nm	% RSD < 1.0	0.27	%
	Mn	257.610 nm	% RSD < 1.0	0.41	%
	La	379.478 nm	% RSD < 1.0	0.57	%
	Ba	455.403 nm	% RSD < 1.0	0.33	%
	Ba	493.408 nm	% RSD < 1.0	0.26	%
Detection Limits : Axial	Tl	190.080 nm	3(sd)	5.51	ppb
	As	193.696 nm	3(sd)	8.59	ppb
	Pb	220.353 nm	3(sd)	0.50	ppb
Detection Limits : Radial	As	193.696 nm	3(sd)	21.00	ppb
	Zn	213.856 nm	3(sd)	0.32	ppb
	Mn	257.610 nm	3(sd)	0.18	ppb
	La	379.478 nm	3(sd)	0.44	ppb
	Ba	455.403 nm	3(sd)	0.17	ppb
	Ba	493.408 nm	3(sd)	0.12	ppb
BEC : Axial (IB X 500)/(IS-IB)	Cd	226.502 nm	≤ 150 ppb	12.46	
BEC : Radial (IB X 1000)/(IS-IB)	Mn	257.610 nm	≤ 45 ppb	30.82	



MAINTENANCE AND TEST CERTIFICATE MODEL

OPTIMA 5300DV

SERIAL NUMBER 077C8011701**DATE TESTED** January 12, 2022**Remarks :**

Commissioning follow as commissioning performance sheets.

This is to certify that the above tests have been performed and the configuration tested



meets



does not meet

the PerkinElmer Specifications listed on this certificate.

This certificate does not modify PerkinElmer's standard terms and condition of sale,
including warranty terms.


Service Department PerkinElmer Ltd.

Authorized Represent

Service Engineer

PinAAcle 900Z Preventive Maintenance (PM)

Company Name:	S.P.S. CONSULTING SERVICE CO., LTD		
Address (Instrument Location):	7 Soi Phaholyothin 24 Phaholyothin Road., Jompol, Chatuchak, Bangkok 10110		
Serial Number:	PZAS19090402	PM Number:	2/2
Customer Name (if applicable):	K. PHENPHA	Telephone Number:	083-926-9252
Customer Support Engineer Name:	K. DUANG	Service Order Number:	WO-01473846
Date PM Performed: (DD-MMM-YYYY)	01-Dec-2021	Next PM Due Date: (DD-MMM-YYYY)	01-Jun-2022
Standard Labor Hours to Complete PM :		5 hours	

Part Number	Release	Publication Date	
09370144 Rev.9	A	January 2018	

Scope

The purpose of this PM is to ensure the continued functionality of the PinAAcle 900Z by inspecting and replacing any worn or damaged parts. This service should only be performed by a trained representative of PerkinElmer.

The customer should save their method before the PM begins.

General Instructions:

The customer must provide the engineer operational data to demonstrate recent instrument performance prior to starting the PM.

Always check with the customer before making any changes that may affect the customer's analysis or calibration, including a current back-up of system software and/or data files.

The completed document should be signed by an authorized PerkinElmer and customer representative and left with the customer.

Update the PM sticker and instrument logbook as required.

Copyright Information

This document contains proprietary information that is protected by copyright. All rights are reserved.

No part of this publication may be reproduced in any form whatsoever or translated into any language without the prior, written permission of PerkinElmer, Inc. **Copyright © 2013 PerkinElmer, Inc.**

Trademarks

Registered names, trademarks, etc. used in this document, even when not specifically marked as such, are protected by law. PerkinElmer is a registered trademark of PerkinElmer, Inc. All other trademarks and registered trademarks not owned by PerkinElmer, Inc. or its subsidiaries that are depicted herein are the property of their respective owners. **Except as specifically set forth in its terms and conditions of sale, PerkinElmer makes no Warranty of any kind with regard to this document, including, but not limited to, the implied warranties of merchantability and fitness for a particular purpose.**

PerkinElmer shall not be liable for incidental or consequential damages in connection with the furnishing or use of this document.

Component List

Component / Specific Model	Serial #	Configuration Notes

Parts Lists

Parts Included with the PM		
Part Number (if applicable)	Description	Quantity
B0501696	Fan Filters	N/A
B3002013	THGA Contact Cylinders	N/A
B3141064	Glycerol for THGA Cooling	N/A

Additional Reagents and Standards Required for PM				
Part Number (if applicable)	Description	Quality	Batch/Lot #	Expired Date (MM/YY)
N9300244	GFAAS Mixed Standard	AR	53-255CRY1	28-Feb-2022

Additional Reagents and Standards Required for PM (Customer Support Solution)				
Part Number (if applicable)	Description	Quantity	Batch/Lot #	Expiration Date (MM/YY)
N/A	DI Water	250 ml.	AR	AR
N/A	0.5% HNO ₃	250 ml.	AR	AR

Additional Tools Required for PM			
Part Number (if applicable)	Description	Quantity	Serial #
B3100652 Or N9307029	Electronic Flow Meter	1	NA
B0505495	Test Jig	1	NA
03030997	System 2 EDL Driver	1	03030997
N3050605	As System 2 EDL	1	16148
N3050121	Cu Lumina HCL	1	092216-010130
N3050109	Ba Lumina HCL	1	102416-040160
N3050139	K Lumina HCL	1	110716-010060
N3050152	Ni Lumina HCL	1	100516-030190
N3050119	Cr Lumina HCL	1	091911-020150

Procedure Checklist

Use (✓) to check off those steps in the checklist that have been completed.

1. General:

- ☒ Review the instrument performance with the customer and document any recent problems.
- ☒ Inspect the customer log book and make any appropriate PM entries.
- ☒ Perform general inspection of system for cleanliness.

2. PC Instrument Software:

- ☒ Instrument Software user files/databases archived, packed, and/or deleted as needed.

3. Mechanical:

- ☒ Inspect and clean all fans and filters. Replace filters if necessary
- ☒ Inspect all gas and water lines for leaks and/or wear. Replace if needed. Thoroughly inspect all quick connects. Replace the Y connector, P/N 09921079, if needed.
- ☒ Clean exterior of the instrument.
- ☒ Check the drain system for signs of wear. Replace worn or damaged parts.
- ☒ Inspect the pole pieces and clean where the pole pieces contact the furnace. Replace the pole piece p-rings as needed, P/N's B0501018 & B0501250. Grease the O-rings as needed with Apiezon L grease, P/N 09905148
- ☒ Inspect the four insulation pads on the front contact housing of the THGA in furnace. If the pads are missing replace the THGA furnace or replace the insulator pads on the furnace.
- ☒ Inspect the graphite tube and clean the contact cylinders. Replace if necessary.
- ☒ Check internal and external gas flows with the Electronic Gas Flow Meter and the Gas Flow Test Probe as described in the Service Manual. Correct if necessary.
- ☒ Check furnace open/close function.
- ☒ Verify the operation of the GFTV Camera for proper operation and viewing alignment in the furnace camera Tube View window. Align if needed.
- ☒ Check the operation of the Halogen Light ASSY for the GFTV Camera. Replace if needed.
- ☒ Check the water level/quality in the recirculation (if applicable). Add distilled water if necessary.
- ☒ Check the cooling system fluid flow rate with the FCS In-Line Flow Meter for proper levels if needed. Refer to SDB# COSY008.STN
- ☒ Perform Cooling System maintenance if needed per SDB# COSY005.STN.
- ☒ Check auto sampler operation.
- ☒ Perform an auto sampler check valve test as described in the Service Manual.
- ☒ Lubricate the spindles of the auto sampler pumps and all moving parts of the tray mechanics as described in the Service Manual.
- ☒ Inspect the auto sampler sampling capillary as described in the Service Manual. Replace if necessary.
- ☒ Inspect the four insulation pads on the front contact housing of the THGA in furnace. If the pads are missing replace the THGA furnace or replace the insulator pads on the furnace.
- ☒ Inspect the graphite tube and clean the contact cylinders. Replace if necessary.
- ☒ Check internal and external gas flows with the Electronic Gas Flow Meter and the Gas Flow Test Probe as described in the Service Manual. Correct if necessary.
- ☒ Check furnace open/close function

4. Electrical:

- ✓ Inspect PC boards. Clean if necessary.
- ✓ Check instrument firmware revisions upgrade to current levels (if necessary)
- ✓ Run Diagnostics Test within the Advanced function of the Spectrometer page. Check the results in the service log folder in the Spectrometer BM Log Viewer.

5. Optics:

- ✓ Inspect and clean the sample compartment windows, if needed.
- ✓ Inspect and clean the furnace windows, if needed.
- ✓ Inspect and clean the GFTV camera lens, if needed.
- ✓ Inspect optics. Clean or replace if necessary,

6. Gasses:

- ✓ Verify that the Gasses supplied to the instrument are within the pressure and purity specifications found in the PinAAcle 900 Series Pre-installation Checklist SDB.
- ✓ Verify that the air filter element is dry. Replace if necessary.

7. After PM Performance tests [THGA]:

7.1 Furnace Gas Flows

Description: Ensures the flow rates are within specification.

Parameter	Specification	Test Results	Pass/Fail
Internal Flow Rate	250 mL/min \pm 25 mL/min	250	Passed
External Flow Rate	100 mL/min \pm 10 mL/min	100	Passed

7.2 Chromium Baseline Noise

Description: Signal to noise check.

Parameter	Specification	Results	Pass/Fail
Baseline Noise	\leq 0.005 Abs.	0.0001	Passed
Standard Deviation	\leq 0.005	0.0001	Passed

7.3 Chromium Characteristic Mass and Precision

Description: Calculate the characteristic mass using the characteristic mass tool and precision from the integrated absorbance values.

Parameter	Specification	Results	Pass/Fail
Cr m ₀ Results	\leq 7.0 pg/0.0044 A-s	3.8	Passed
Precision	\leq 2.0 %	1.64	Passed

7.4 Copper Characteristic Mass and Zeeman Ratio

Description: Calculate the characteristic mass using the characteristic mass tool and check the Zeeman Ratio.

Parameter	Specification	Results	Pass/Fail
Cu m ₀ Result	≤ 16.5 pg/0.0044 A-s	13.9	Passed
Zeeman Ratio	0.52 ± 0.04	0.52	Passed

8. Review:

- ☒ Review with the customer PM work performed.
- ☒ Review with the customer routine maintenance procedures.
- ☒ Discuss recommended customer supplied materials to have on hand.
- ☒ Attach PM sticker.

Additional Comments

Additional Comments Regarding the PM	
Zeeman Ratio	$= \frac{\text{Atomic Signal (Peak area)}}{\text{Atomic Signal (Peak area)} + \text{Background Signal (Peak area)}}$ $= \frac{0.1593}{0.1593+0.1414}$ $= 0.52$

Review

<p><i>The preventive maintenance checks and if applicable performance tests for PinAAcle 900Z have been completed.</i></p> <p><i>This PinAAcle 900Z Passes <input checked="" type="checkbox"/> Fails <input type="checkbox"/> the preventive maintenance.</i></p>		
Review of Preventive Maintenance:		
Authorized PerkinElmer Representative:	<div style="background-color: black; width: 100px; height: 100px;"></div>	Date: 01-Dec-2021 <small>(DD-MMM-YYYY)</small>
Authorized Customer Representative:	<div style="background-color: black; width: 100px; height: 100px;"></div>	Date: 01-Dec-2021 <small>(DD-MMM-YYYY)</small>