

คุณภาพน้ำผิวดิน



ANALYSIS REPORT

CUSTOMER NAME : ECO ORIENT RESOURCES (THAILAND) LTD.
ADDRESS : 555 RASA TOWER II, 12TH FLOOR, UNIT 1203 PHAHOLYOTHIN ROAD CHATUCHAK CHATUCHAK BANGKOK 10900
CONTACT INFORMATION : TEL : 0 2937 1124-9 e-mail : anucha@ecothai.net
SAMPLING SOURCE : SW18: HUI WANG CHAK (UPSTREAM) (UTM WGS 84 ZONE 47P 737104E 1725414N)
SAMPLE TYPE : SURFACE WATER
SAMPLING DATE : FEBRUARY 9, 2022
SAMPLING TIME : 14:20 HOUR
SAMPLING METHOD ° : GRAB, GRAB AND STERILE TECHNIQUE
SAMPLING BY ° : MR TEERAPONG SRIKAMHANG
ANALYZED BY : MISS PORNPIMOL WAENTHONG

RECEIVED DATE : FEBRUARY 10, 2022
ANALYTICAL DATE : FEBRUARY 10-19, 2022
REPORT NO. : 2022-U012753
WORK NO. : 2022-000219
ANALYSIS NO. : T22AC416-0001

PARAMETER	UNIT	METHOD OF ANALYSIS	RESULT	REGULATORY STANDARD	DETECTION LIMIT
			SURFACE WATER T22AC416-0001		
pH °	-	ELECTROMETRIC METHOD AT SITE (SM:4500-H° B)	7.5 (30°C)	5.0-9.0	-
ELECTRICAL CONDUCTIVITY °	µmhos/cm	ELECTRICAL CONDUCTIVITY METHOD AT SITE (SM: 2510 B)	143 (30°C)	-	0.1
SULPHATE °	mg/L SO ₄ ²⁻	TURBIDIMETRIC METHOD (SM: 4500-SO ₄ ²⁻ E)	2.2	-	0.3
TOTAL PETROLEUM HYDROCARBONS °	mg/L	SOXHLET EXTRACTION METHOD (SM: 5520 D AND 5520 F)	ND	-	3
METALS					
ARSENIC °	mg/L As	HYDRIDE GENERATION AAS METHOD (SM: 3114 C)	ND	≤ 0.01	0.0003
CADMIUM °	mg/L Cd	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.005*, ≤ 0.05**	0.002
CHROMIUM °	mg/L Cr	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	-	0.005
COPPER °	mg/L Cu	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	< LOQ	≤ 0.1	0.002
IRON °	mg/L Fe	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	2.98	-	0.005
LEAD °	mg/L Pb	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.05	0.003
MANGANESE °	mg/L Mn	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	0.280	≤ 1.0	0.002
MERCURY °	mg/L Hg	IN-HOUSE METHOD: UAE.TP.HEM.002 (COLD VAPOUR ATOMIC ABSORPTION SPECTROMETRIC METHOD); SM: 3112 B	ND	≤ 0.002	0.0001
NICKEL °	mg/L Ni	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.1	0.005



PARAMETER	UNIT	METHOD OF ANALYSIS	RESULT	REGULATORY STANDARD	DETECTION LIMIT
			SURFACE WATER T22AC416-0001		
ZINC ^c	mg/L Zn	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 1.0	0.003
SAMPLE CONDITION WATER'S COLOUR/TURBID SEDIMENT			YELLOW/TURBID BROWN		

^a : ISO/IEC 17025 ACCREDITED BY THAI INDUSTRIAL STANDARDS INSTITUTE (TISI)

^b : ISO/IEC 17025 ACCREDITED BY DEPARTMENT OF SCIENCE SERVICE (DSS)

^c : VERIFIED BY OWN LABORATORY QUALITY SYSTEM, BUT STILL NOT ACCREDITED

IN-HOUSE : BASED ON STANDARD METHODS FOR THE EXAMINATION OF WATER AND WASTEWATER, APHA, AWWA, WEF, 23rd EDITION, 2017.

SM : STANDARD METHODS FOR THE EXAMINATION OF WATER AND WASTEWATER, APHA, AWWA, WEF, 23rd EDITION, 2017.

REGULATORY STANDARD : SURFACE WATER QUALITY STANDARDS CLASS 3, NOTIFICATION OF THE NATIONAL ENVIRONMENT BOARD, NO.8, B.E. 2537 ISSUED UNDER THE ENCHANCEMENT AND CONSERVATION OF NATIONAL ENVIRONMENTAL QUALITY ACT, B.E. 2535, PUBLISHED IN THE ROYAL GOVERNMENT GAZETTE, VOL. 111, PART 16, DATED FEBRUARY 24, B.E. 2537 (1994).

CLASS 3 : MEDIUM CLEAN FRESH SURFACE WATER RESOURCES USED FOR

- (1) CONSUMPTION, BUT PASSING THROUGH ON ORDINARY TREATMENT PROCESS BEFORE USING
- (2) AGRICULTURE

≤ 0.005* : WHEN WATER HARDNESS NOT MORE THAN 100 mg/L AS CaCO₃

≤ 0.05** : WHEN WATER HARDNESS MORE THAN 100 mg/L AS CaCO₃

ND : NON-DETECTABLE.

< LOQ : < LEVEL OF QUANTITATION (COPPER ≥ 0.002 AND < 0.025 mg/L).

Benjawan V.

(MISS BENJAWAN VIRIYOTHAJ)
LABORATORY SUPERVISOR

FEBRUARY 28, 2022

ANALYSIS REPORT

CUSTOMER NAME : ECO ORIENT RESOURCES (THAILAND) LTD.
ADDRESS : 555 RASA TOWER II, 12TH FLOOR, UNIT 1203 PHAHOLYOTHIN ROAD CHATUCHAK CHATUCHAK BANGKOK 10900
CONTACT INFORMATION : TEL : 0 2937 1124-9 e-mail : anucha@ecothai.net
SAMPLING SOURCE : SW19: HUI WANG CHAK (DOWNSTREAM) (UTM WGS 84 ZONE 47P 735421E 1726710N)
SAMPLE TYPE : SURFACE WATER
SAMPLING DATE : FEBRUARY 9, 2022
SAMPLING TIME : 14:55 HOUR
SAMPLING METHOD ° : GRAB, GRAB AND STERILE TECHNIQUE
SAMPLING BY ° : MR TEERAPONG SRIKAMHANG
ANALYZED BY : MISS PORNPIMOL WAENTHONG

RECEIVED DATE : FEBRUARY 10, 2022
ANALYTICAL DATE : FEBRUARY 10-19, 2022
REPORT NO. : 2022-U012754
WORK NO. : 2022-000219
ANALYSIS NO. : T22AC416-0002

PARAMETER	UNIT	METHOD OF ANALYSIS	RESULT	REGULATORY STANDARD	DETECTION LIMIT
			SURFACE WATER T22AC416-0002		
pH °	-	ELECTROMETRIC METHOD AT SITE (SM:4500-H ⁺ B)	7.5 (30°C)	5.0-9.0	-
ELECTRICAL CONDUCTIVITY °	µmhos/cm	ELECTRICAL CONDUCTIVITY METHOD AT SITE (SM: 2510 B)	124 (30°C)	-	0.1
SULPHATE °	mg/L SO ₄ ²⁻	TURBIDIMETRIC METHOD (SM: 4500-SO ₄ ²⁻ E)	ND	-	0.3
TOTAL PETROLEUM HYDROCARBONS °	mg/L	SOXHLET EXTRACTION METHOD (SM: 5520 D AND 5520 F)	ND	-	3
METALS					
ARSENIC °	mg/L As	HYDRIDE GENERATION AAS METHOD (SM: 3114 C)	0.0003	≤ 0.01	0.0003
CADMIUM °	mg/L Cd	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.005*, ≤ 0.05**	0.002
CHROMIUM °	mg/L Cr	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	-	0.005
COPPER °	mg/L Cu	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.1	0.002
IRON °	mg/L Fe	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	2.04	-	0.005
LEAD °	mg/L Pb	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.05	0.003
MANGANESE °	mg/L Mn	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	1.05	≤ 1.0	0.002
MERCURY ^b	mg/L Hg	IN-HOUSE METHOD: UAE.TP.HEM.002 (COLD VAPOUR ATOMIC ABSORPTION SPECTROMETRIC METHOD); SM: 3112 B	ND	≤ 0.002	0.0001
NICKEL °	mg/L Ni	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.1	0.005



PARAMETER	UNIT	METHOD OF ANALYSIS	RESULT	REGULATORY STANDARD	DETECTION LIMIT
			SURFACE WATER T22AC416-0002		
ZINC ^c	mg/L Zn	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 1.0	0.003
SAMPLE CONDITION WATER'S COLOUR/TURBID SEDIMENT			YELLOW/CLEAR BROWN		

^a : ISO/IEC 17025 ACCREDITED BY THAI INDUSTRIAL STANDARDS INSTITUTE (TISI)

^b : ISO/IEC 17025 ACCREDITED BY DEPARTMENT OF SCIENCE SERVICE (DSS)

^c : VERIFIED BY OWN LABORATORY QUALITY SYSTEM, BUT STILL NOT ACCREDITED

IN-HOUSE : BASED ON STANDARD METHODS FOR THE EXAMINATION OF WATER AND WASTEWATER, APHA, AWWA, WEF, 23rd EDITION, 2017.

SM : STANDARD METHODS FOR THE EXAMINATION OF WATER AND WASTEWATER, APHA, AWWA, WEF, 23rd EDITION, 2017.

REGULATORY STANDARD : SURFACE WATER QUALITY STANDARDS CLASS 3, NOTIFICATION OF THE NATIONAL ENVIRONMENT BOARD, NO.8, B.E. 2537 ISSUED UNDER THE ENCHANCEMENT AND CONSERVATION OF NATIONAL ENVIRONMENTAL QUALITY ACT, B.E. 2535, PUBLISHED IN THE ROYAL GOVERNMENT GAZETTE, VOL. 111, PART 16, DATED FEBRUARY 24, B.E. 2537 (1994).

CLASS 3 : MEDIUM CLEAN FRESH SURFACE WATER RESOURCES USED FOR
(1) CONSUMPTION, BUT PASSING THROUGH ON ORDINARY TREATMENT PROCESS BEFORE USING
(2) AGRICULTURE

≤ 0.005* : WHEN WATER HARDNESS NOT MORE THAN 100 mg/L AS CaCO₃

≤ 0.05** : WHEN WATER HARDNESS MORE THAN 100 mg/L AS CaCO₃

ND : NON-DETECTABLE.

Benjawan V.

(MISS BENJAWAN VIRIYOTHAI)
LABORATORY SUPERVISOR

FEBRUARY 28, 2022

ระยะผลิตปิโตรเลียม



คุณภาพน้ำใต้ดิน



ANALYSIS REPORT

CUSTOMER NAME : ECO ORIENT RESOURCES (THAILAND) LTD.
ADDRESS : 555 RASA TOWER II, 12TH FLOOR, UNIT 1203 PHAHOLYOTHIN ROAD CHATUCHAK CHATUCHAK BANGKOK 10900
CONTACT INFORMATION : TEL : 0 2937 1124-9 e-mail : anucha@eco-thai.net
SAMPLING SOURCE : MWBRN-5 (UP GRADIENT) (UTM WGS 84 ZONE 47P 736142E 1728147N)
SAMPLE TYPE : GROUNDWATER
SAMPLING DATE : MARCH 30, 2022
SAMPLING TIME : 15:30 HOUR
SAMPLING METHOD ° : SUBMERSIBLE PUMP
SAMPLING BY ° : MR PORAWORN BUNNAG
ANALYZED BY : MISS NADNAPA KAMOLBOON

RECEIVED DATE : APRIL 1, 2022
ANALYTICAL DATE : APRIL 1-20, 2022
REPORT NO. : 2022-U028862
WORK NO. : 2022-001206
ANALYSIS NO. : T22AG205-0001

PARAMETER	UNIT	METHOD OF ANALYSIS	RESULT	REGULATORY STANDARD	DETECTION LIMIT
			MWBRN-5 (UP GRADIENT) T22AG205-0001		
pH °	-	ELECTROMETRIC METHOD AT SITE (SM:4500-H+ B)	8.0 (30°C)	-	-
CONDUCTIVITY °	µmhos/cm	ELECTRICAL CONDUCTIVITY METHOD AT SITE (SM: 2510 B)	997 (30°C)	-	0.1
SULPHATE °	mg/L SO ₄ ²⁻	TURBIDIMETRIC METHOD (SM: 4500-SO ₄ ²⁻ E)	ND	-	0.3
TOTAL PETROLEUM HYDROCARBONS °	mg/L	SOXHLET EXTRACTION METHOD (SM: 5520 D AND 5520 F)	ND	-	3
METALS					
ARSENIC °	mg/L As	HYDRIDE GENERATION AAS METHOD (SM: 3114 C)	0.0015	≤ 0.01	0.0003
IRON °	mg/L Fe	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	0.236	-	0.005
TOTAL CHROMIUM °	mg/L Cr	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	-	0.005
MERCURY °	mg/L Hg	IN-HOUSE METHOD: UAE.TP.HEM.002 (COLD VAPOUR ATOMIC ABSORPTION SPECTROMETRIC METHOD); SM: 3112 B	ND	≤ 0.001	0.0001
CADMIUM °	mg/L Cd	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.003	0.002
COPPER °	mg/L Cu	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 1.0	0.002
LEAD °	mg/L Pb	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.01	0.003
MANGANESE °	mg/L Mn	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	0.255	≤ 0.5	0.002



PARAMETER	UNIT	METHOD OF ANALYSIS	RESULT	REGULATORY STANDARD	DETECTION LIMIT
			MWBRN-5 (UP GRADIENT) T22AG205-0001		
NICKEL ^c	mg/L Ni	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.02	0.005
ZINC ^c	mg/L Zn	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 5.0	0.003
SAMPLE CONDITION WATER'S COLOUR/TURBID SEDIMENT			COLOURLESS/CLEAR		

^a : ISO/IEC 17025 ACCREDITED BY THAI INDUSTRIAL STANDARDS INSTITUTE (TISI)

^b : ISO/IEC 17025 ACCREDITED BY DEPARTMENT OF SCIENCE SERVICE (DSS)

^c : VERIFIED BY OWN LABORATORY QUALITY SYSTEM, BUT STILL NOT ACCREDITED

IN-HOUSE : BASED ON STANDARD METHODS FOR THE EXAMINATION OF WATER AND WASTEWATER, APHA, AWWA, WEF, 23rd EDITION, 2017.

SM : STANDARD METHODS FOR THE EXAMINATION OF WATER AND WASTEWATER, APHA, AWWA, WEF, 23rd EDITION, 2017.

REGULATORY STANDARD : GROUNDWATER QUALITY STANDARDS, NOTIFICATION OF THE NATIONAL ENVIRONMENT BOARD NO.20 (B.E. 2543)
ISSUED UNDER THE ENHANCEMENT AND CONSERVATION OF NATIONAL ENVIRONMENTAL QUALITY ACT B.E. 2535.

ND : NON-DETECTABLE.

Benjawan V.

(MISS BENJAWAN VIRIYOTHAI)
LABORATORY SUPERVISOR

APRIL 22, 2022

ANALYSIS REPORT

CUSTOMER NAME : ECO ORIENT RESOURCES (THAILAND) LTD.
ADDRESS : 555 RASA TOWER II, 12TH FLOOR, UNIT 1203 PHAHOLYOTHIN ROAD CHATUCHAK CHATUCHAK BANGKOK 10900
CONTACT INFORMATION : TEL : 0 2937 1124-9 e-mail : anucha@eco-thai.net
SAMPLING SOURCE : MWBRN-5 (DOWN GRADIENT) (UTM WGS 84 ZONE 47P 736000E 1728040N)
SAMPLE TYPE : GROUNDWATER
SAMPLING DATE : MARCH 30, 2022
SAMPLING TIME : 16:00 HOUR
SAMPLING METHOD : SUBMERSIBLE PUMP
SAMPLING BY : MR PORAWORN BUNNAG
ANALYZED BY : MISS NADNAPA KAMOLBOON

RECEIVED DATE : APRIL 1, 2022
ANALYTICAL DATE : APRIL 1-20, 2022
REPORT NO. : 2022-U028863
WORK NO. : 2022-001206
ANALYSIS NO. : T22AG205-0002

PARAMETER	UNIT	METHOD OF ANALYSIS	RESULT	REGULATORY STANDARD	DETECTION LIMIT
			MWBRN-5 (DOWN GRADIENT) T22AG205-0002		
pH ^c	-	ELECTROMETRIC METHOD AT SITE (SM:4500-H ⁺ B)	8.5 (28°C)	-	-
CONDUCTIVITY ^c	µmhos/cm	ELECTRICAL CONDUCTIVITY METHOD AT SITE (SM: 2510 B)	413 (28°C)	-	0.1
SULPHATE ^c	mg/L SO ₄ ²⁻	TURBIDIMETRIC METHOD (SM: 4500-SO ₄ ²⁻ E)	0.9	-	0.3
TOTAL PETROLEUM HYDROCARBONS ^c	mg/L	SOXHLET EXTRACTION METHOD (SM: 5520 D AND 5520 F)	ND	-	3
METALS					
ARSENIC ^c	mg/L As	HYDRIDE GENERATION AAS METHOD (SM: 3114 C)	0.0005	≤ 0.01	0.0003
IRON ^c	mg/L Fe	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	0.479	-	0.005
TOTAL CHROMIUM ^c	mg/L Cr	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	-	0.005
MERCURY ^b	mg/L Hg	IN-HOUSE METHOD: UAE.TP.HEM.002 (COLD VAPOUR ATOMIC ABSORPTION SPECTROMETRIC METHOD); SM: 3112 B	ND	≤ 0.001	0.0001
CADMIUM ^c	mg/L Cd	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.003	0.002
COPPER ^c	mg/L Cu	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 1.0	0.002
LEAD ^c	mg/L Pb	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	0.136	≤ 0.01	0.003
MANGANESE ^c	mg/L Mn	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	0.434	≤ 0.5	0.002



PARAMETER	UNIT	METHOD OF ANALYSIS	RESULT	REGULATORY STANDARD	DETECTION LIMIT
			MWBRN-5 (DOWN GRADIENT) T22AG205-0002		
NICKEL ^c	mg/L Ni	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.02	0.005
ZINC ^c	mg/L Zn	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 5.0	0.003
SAMPLE CONDITION WATER'S COLOUR/TURBID SEDIMENT			COLOURLESS/CLEAR		

^a : ISO/IEC 17025 ACCREDITED BY THAI INDUSTRIAL STANDARDS INSTITUTE (TISI)

^b : ISO/IEC 17025 ACCREDITED BY DEPARTMENT OF SCIENCE SERVICE (DSS)

^c : VERIFIED BY OWN LABORATORY QUALITY SYSTEM, BUT STILL NOT ACCREDITED

IN-HOUSE : BASED ON STANDARD METHODS FOR THE EXAMINATION OF WATER AND WASTEWATER, APHA, AWWA, WEF, 23rd EDITION, 2017.

SM : STANDARD METHODS FOR THE EXAMINATION OF WATER AND WASTEWATER, APHA, AWWA, WEF, 23rd EDITION, 2017.

REGULATORY STANDARD : GROUNDWATER QUALITY STANDARDS, NOTIFICATION OF THE NATIONAL ENVIRONMENT BOARD NO.20 (B.E. 2543)
ISSUED UNDER THE ENHANCEMENT AND CONSERVATION OF NATIONAL ENVIRONMENTAL QUALITY ACT B.E. 2535.

ND : NON-DETECTABLE.

Benjawan V.

(MISS BENJAWAN VIRIYOTHAI)
LABORATORY SUPERVISOR

APRIL 22, 2022

ANALYSIS REPORT

CUSTOMER NAME : ECO ORIENT RESOURCES (THAILAND) LTD.
ADDRESS : 555 RASA TOWER II, 12TH FLOOR, UNIT 1203 PHAHOLYOTHIN ROAD CHATUCHAK CHATUCHAK BANGKOK 10900
CONTACT INFORMATION : TEL : 0 2937 1124-9 e-mail : anucha@ecothai.net
SAMPLING SOURCE : BRN-5-GW1 (UTM WGS 84 ZONE 47P 737021E 1729133N)
SAMPLE TYPE : GROUNDWATER
SAMPLING DATE : MARCH 29, 2022
SAMPLING TIME : 15:55 HOUR
SAMPLING METHOD : GRAB
SAMPLING BY : MR PORAWORN BUNNAG ๓-145-๓-0102
ANALYZED BY : MISS CHOMTHANAN APHIPATPAPHA ๓-145-๓-0022

RECEIVED DATE : MARCH 30, 2022
ANALYTICAL DATE : MARCH 30-APRIL 14, 2022
REPORT NO. : 2022-U027995
WORK NO. : 2022-001206
ANALYSIS NO. : T22AG078-0001

PARAMETER	UNIT	METHOD OF ANALYSIS	RESULT	REGULATORY STANDARD	DETECTION LIMIT
			BRN-5-GW1 T22AG078-0001		
pH ^c	-	ELECTROMETRIC METHOD AT SITE (SM:4500-H ⁺ B)	7.6 (29°C)	-	-
CONDUCTIVITY ^c	µmhos/cm	ELECTRICAL CONDUCTIVITY METHOD AT SITE (SM: 2510 B)	201 (29°C)	-	0.1
SULPHATE ^c	mg/L SO ₄ ²⁻	TURBIDIMETRIC METHOD (SM: 4500-SO ₄ ²⁻ E)	4.5	-	0.3
TOTAL PETROLEUM HYDROCARBONS ^c	mg/L	SOXHLET EXTRACTION METHOD (SM: 5520 D AND 5520 F)	ND	-	3
METALS					
ARSENIC ^c	mg/L As	HYDRIDE GENERATION AAS METHOD (SM: 3114 C)	ND	≤ 0.01	0.0003
IRON ^c	mg/L Fe	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	0.269	-	0.005
TOTAL CHROMIUM ^c	mg/L Cr	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	-	0.005
MERCURY ^b	mg/L Hg	IN-HOUSE METHOD: UAE.TP.HEM.002 (COLD VAPOUR ATOMIC ABSORPTION SPECTROMETRIC METHOD); SM: 3112 B	< LOQ	≤ 0.001	0.0001
CADMIUM ^c	mg/L Cd	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.003	0.002
COPPER ^c	mg/L Cu	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 1.0	0.002
LEAD ^c	mg/L Pb	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.01	0.003
MANGANESE ^c	mg/L Mn	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	0.052	≤ 0.5	0.002



PARAMETER	UNIT	METHOD OF ANALYSIS	RESULT	REGULATORY STANDARD	DETECTION LIMIT
			BRN-5-GW1 T22AG078-0001		
NICKEL ^c	mg/L Ni	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.02	0.005
ZINC ^c	mg/L Zn	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	< LOQ	≤ 5.0	0.003
SAMPLE CONDITION WATER'S COLOUR/TURBID SEDIMENT			COLOURLESS/CLEAR		

^a : ISO/IEC 17025 ACCREDITED BY THAI INDUSTRIAL STANDARDS INSTITUTE (TISI)

^b : ISO/IEC 17025 ACCREDITED BY DEPARTMENT OF SCIENCE SERVICE (DSS)

^c : VERIFIED BY OWN LABORATORY QUALITY SYSTEM, BUT STILL NOT ACCREDITED

IN-HOUSE : BASED ON STANDARD METHODS FOR THE EXAMINATION OF WATER AND WASTEWATER, APHA, AWWA, WEF, 23rd EDITION, 2017.

SM : STANDARD METHODS FOR THE EXAMINATION OF WATER AND WASTEWATER, APHA, AWWA, WEF, 23rd EDITION, 2017.

REGULATORY STANDARD : GROUNDWATER QUALITY STANDARDS, NOTIFICATION OF THE NATIONAL ENVIRONMENT BOARD NO.20 (B.E. 2543)
ISSUED UNDER THE ENHANCEMENT AND CONSERVATION OF NATIONAL ENVIRONMENTAL QUALITY ACT B.E. 2535.

ND : NON-DETECTABLE.

< LOQ : < LEVEL OF QUANTITATION (MERCURY ≥ 0.0001 AND < 0.0005 mg/L, ZINC ≥ 0.003 AND < 0.025 mg/L).

Benjawan V.

(MISS BENJAWAN VIRIYOTHAI)

LABORATORY SUPERVISOR

๓-145-๓-0006

APRIL 20, 2022

ANALYSIS REPORT

CUSTOMER NAME : ECO ORIENT RESOURCES (THAILAND) LTD.
ADDRESS : 555 RASA TOWER II, 12TH FLOOR, UNIT 1203 PHAHOLYOTHIN ROAD CHATUCHAK CHATUCHAK BANGKOK 10900
CONTACT INFORMATION : TEL : 0 2937 1124-9 e-mail : anucha@ecothai.net
SAMPLING SOURCE : BRN-5-GW2 (UTM WGS 84 ZONE 47P 734949E 1727784N)
SAMPLE TYPE : GROUNDWATER
SAMPLING DATE : MARCH 29, 2022
SAMPLING TIME : 15:40 HOUR
SAMPLING METHOD : GRAB
SAMPLING BY : MR PORAWORN BUNNAG
ANALYZED BY : MISS CHOMTHANAN APHIPATPAPHA

RECEIVED DATE : MARCH 30, 2022
ANALYTICAL DATE : MARCH 30-APRIL 14, 2022
REPORT NO. : 2022-U027996
WORK NO. : 2022-001206
ANALYSIS NO. : T22AG078-0002

PARAMETER	UNIT	METHOD OF ANALYSIS	RESULT	REGULATORY STANDARD	DETECTION LIMIT
			BRN-5-GW2 T22AG078-0002		
pH ^c	-	ELECTROMETRIC METHOD AT SITE (SM:4500-H ⁺ B)	7.1 (32°C)	-	-
CONDUCTIVITY ^c	µmhos/cm	ELECTRICAL CONDUCTIVITY METHOD AT SITE (SM: 2510 B)	277 (32°C)	-	0.1
SULPHATE ^c	mg/L SO ₄ ²⁻	TURBIDIMETRIC METHOD (SM: 4500-SO ₄ ²⁻ E)	ND	-	0.3
TOTAL PETROLEUM HYDROCARBONS ^c	mg/L	SOXHLET EXTRACTION METHOD (SM: 5520 D AND 5520 F)	ND	-	3
METALS					
ARSENIC ^c	mg/L As	HYDRIDE GENERATION AAS METHOD (SM: 3114 C)	ND	≤ 0.01	0.0003
IRON ^c	mg/L Fe	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	0.169	-	0.005
TOTAL CHROMIUM ^c	mg/L Cr	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	-	0.005
MERCURY ^b	mg/L Hg	IN-HOUSE METHOD: UAE.TP.HEM.002 (COLD VAPOUR ATOMIC ABSORPTION SPECTROMETRIC METHOD); SM: 3112 B	< LOQ	≤ 0.001	0.0001
CADMIUM ^c	mg/L Cd	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.003	0.002
COPPER ^c	mg/L Cu	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	< LOQ	≤ 1.0	0.002
LEAD ^c	mg/L Pb	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.01	0.003
MANGANESE ^c	mg/L Mn	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	0.127	≤ 0.5	0.002



PARAMETER	UNIT	METHOD OF ANALYSIS	RESULT	REGULATORY STANDARD	DETECTION LIMIT
			BRN-5-GW2 T22AG078-0002		
NICKEL ^c	mg/L Ni	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.02	0.005
ZINC ^c	mg/L Zn	IN-HOUSE METHOD: UAE.TP.GW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	< LOQ	≤ 5.0	0.003
SAMPLE CONDITION WATER'S COLOUR/TURBID SEDIMENT			COLOURLESS/CLEAR		

^a : ISO/IEC 17025 ACCREDITED BY THAI INDUSTRIAL STANDARDS INSTITUTE (TISI)

^b : ISO/IEC 17025 ACCREDITED BY DEPARTMENT OF SCIENCE SERVICE (DSS)

^c : VERIFIED BY OWN LABORATORY QUALITY SYSTEM, BUT STILL NOT ACCREDITED

IN-HOUSE : BASED ON STANDARD METHODS FOR THE EXAMINATION OF WATER AND WASTEWATER, APHA, AWWA, WEF, 23rd EDITION, 2017.

SM : STANDARD METHODS FOR THE EXAMINATION OF WATER AND WASTEWATER, APHA, AWWA, WEF, 23rd EDITION, 2017.

REGULATORY STANDARD : GROUNDWATER QUALITY STANDARDS, NOTIFICATION OF THE NATIONAL ENVIRONMENT BOARD NO.20 (B.E. 2543)
ISSUED UNDER THE ENHANCEMENT AND CONSERVATION OF NATIONAL ENVIRONMENTAL QUALITY ACT B.E. 2535.

ND : NON-DETECTABLE.

< LOQ : < LEVEL OF QUANTITATION (MERCURY ≥ 0.0001 AND < 0.0005 mg/L, COPPER ≥ 0.002 AND < 0.025 mg/L,
ZINC ≥ 0.003 AND < 0.025 mg/L).

Benjawan V.

(MISS BENJAWAN VIRIYOTHAJ)
LABORATORY SUPERVISOR

APRIL 20, 2022

คุณภาพน้ำผิวดิน



ANALYSIS REPORT

CUSTOMER NAME : ECO ORIENT RESOURCES (THAILAND) LTD.
ADDRESS : 555 RASA TOWER II, 12TH FLOOR, UNIT 1203 PHAHOLYOTHIN ROAD CHATUCHAK CHATUCHAK BANGKOK 10900
CONTACT INFORMATION : TEL : 0 2937 1124-9 e-mail : anucha@ecothai.net
SAMPLING SOURCE : BRN-5-SW1 (UTM WGS 84 ZONE 47P 736529E 1727573N)
SAMPLE TYPE : SURFACE WATER
SAMPLING DATE : MARCH 31, 2022
SAMPLING TIME : 14:30 HOUR
SAMPLING METHOD ° : GRAB
SAMPLING BY ° : MR PORAWORN BUNNAG
ANALYZED BY : MISS NADNAPA KAMOLBOON

RECEIVED DATE : APRIL 1, 2022
ANALYTICAL DATE : APRIL 1-14, 2022
REPORT NO. : 2022-U028673
WORK NO. : 2022-001206
ANALYSIS NO. : T22AG233-0001

PARAMETER	UNIT	METHOD OF ANALYSIS	RESULT	REGULATORY STANDARD	DETECTION LIMIT
			BRN-5-SW1 T22AG233-0001		
pH °	-	ELECTROMETRIC METHOD AT SITE (SM:4500-H+ B)	7.5 (30°C)	5.0-9.0	-
ELECTRICAL CONDUCTIVITY °	µmhos/cm	ELECTRICAL CONDUCTIVITY METHOD AT SITE (SM: 2510 B)	71.3 (30°C)	-	0.1
SULPHATE °	mg/L SO ₄ ²⁻	TURBIDIMETRIC METHOD (SM: 4500-SO ₄ ²⁻ E)	0.3	-	0.3
TOTAL PETROLEUM HYDROCARBONS °	mg/L	SOXHLET EXTRACTION METHOD (SM: 5520 D AND 5520 F)	ND	-	3
METALS					
ARSENIC °	mg/L As	HYDRIDE GENERATION AAS METHOD (SM: 3114 C)	0.0005	≤ 0.01	0.0003
CADMIUM °	mg/L Cd	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.005*, ≤ 0.05**	0.002
COPPER °	mg/L Cu	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.1	0.002
IRON °	mg/L Fe	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	3.93	-	0.005
LEAD °	mg/L Pb	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.05	0.003
MANGANESE °	mg/L Mn	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	0.326	≤ 1.0	0.002
MERCURY °	mg/L Hg	IN-HOUSE METHOD: UAE.TP.HEM.002 (COLD VAPOUR ATOMIC ABSORPTION SPECTROMETRIC METHOD); SM: 3112 B	ND	≤ 0.002	0.0001
NICKEL °	mg/L Ni	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.1	0.005
TOTAL CHROMIUM °	mg/L Cr	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	-	0.005



PARAMETER	UNIT	METHOD OF ANALYSIS	RESULT	REGULATORY STANDARD	DETECTION LIMIT
			BRN-5-SW1 T22AG233-0001		
ZINC ^c	mg/L Zn	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 1.0	0.003
SAMPLE CONDITION WATER'S COLOUR/TURBID SEDIMENT			YELLOW/TURBID BROWN		

^a : ISO/IEC 17025 ACCREDITED BY THAI INDUSTRIAL STANDARDS INSTITUTE (TISI)

^b : ISO/IEC 17025 ACCREDITED BY DEPARTMENT OF SCIENCE SERVICE (DSS)

^c : VERIFIED BY OWN LABORATORY QUALITY SYSTEM, BUT STILL NOT ACCREDITED

IN-HOUSE : BASED ON STANDARD METHODS FOR THE EXAMINATION OF WATER AND WASTEWATER, APHA, AWWA, WEF, 23rd EDITION, 2017.

SM : STANDARD METHODS FOR THE EXAMINATION OF WATER AND WASTEWATER, APHA, AWWA, WEF, 23rd EDITION, 2017.

REGULATORY STANDARD : SURFACE WATER QUALITY STANDARDS CLASS 3, NOTIFICATION OF THE NATIONAL ENVIRONMENT BOARD, NO.8, B.E. 2537 ISSUED UNDER THE ENHANCEMENT AND CONSERVATION OF NATIONAL ENVIRONMENTAL QUALITY ACT, B.E. 2535, PUBLISHED IN THE ROYAL GOVERNMENT GAZETTE, VOL. 111, PART 16, DATED FEBRUARY 24, B.E. 2537 (1994).

CLASS 3 : MEDIUM CLEAN FRESH SURFACE WATER RESOURCES USED FOR

- (1) CONSUMPTION, BUT PASSING THROUGH ON ORDINARY TREATMENT PROCESS BEFORE USING
- (2) AGRICULTURE

≤ 0.005* : WHEN WATER HARDNESS NOT MORE THAN 100 mg/L AS CaCO₃

≤ 0.05** : WHEN WATER HARDNESS MORE THAN 100 mg/L AS CaCO₃

ND : NON-DETECTABLE

Benjawan V.

(MISS BENJAWAN VIRIYOTHAI)
LABORATORY SUPERVISOR

APRIL 22, 2022

ANALYSIS REPORT

CUSTOMER NAME : ECO ORIENT RESOURCES (THAILAND) LTD.
ADDRESS : 555 RASA TOWER II, 12TH FLOOR, UNIT 1203 PHAHOLYOTHIN ROAD CHATUCHAK CHATUCHAK BANGKOK 10900
CONTACT INFORMATION : TEL : 0 2937 1124-9 e-mail : anucha@ecothai.net
SAMPLING SOURCE : BRN-5-SW2 (UTM WGS 84 ZONE 47P 735637E 1727878N)
SAMPLE TYPE : SURFACE WATER
SAMPLING DATE : MARCH 31, 2022
SAMPLING TIME : 14:00 HOUR
SAMPLING METHOD ° : GRAB
SAMPLING BY ° : MR PORAWORN BUNNAG
ANALYZED BY : MISS NADNAPA KAMOLBOON

RECEIVED DATE : APRIL 1, 2022
ANALYTICAL DATE : APRIL 1-14, 2022
REPORT NO. : 2022-U028674
WORK NO. : 2022-001206
ANALYSIS NO. : T22AG233-0002

PARAMETER	UNIT	METHOD OF ANALYSIS	RESULT	REGULATORY STANDARD	DETECTION LIMIT
			BRN-5-SW2 T22AG233-0002		
pH °	-	ELECTROMETRIC METHOD AT SITE (SM:4500-H ⁺ B)	7.6 (31°C)	5.0-9.0	-
ELECTRICAL CONDUCTIVITY °	µmhos/cm	ELECTRICAL CONDUCTIVITY METHOD AT SITE (SM: 2510 B)	69.2 (31°C)	-	0.1
SULPHATE °	mg/L SO ₄ ²⁻	TURBIDIMETRIC METHOD (SM: 4500-SO ₄ ²⁻ E)	4.6	-	0.3
TOTAL PETROLEUM HYDROCARBONS °	mg/L	SOXHLET EXTRACTION METHOD (SM: 5520 D AND 5520 F)	ND	-	3
METALS					
ARSENIC °	mg/L As	HYDRIDE GENERATION AAS METHOD (SM: 3114 C)	0.0007	≤ 0.01	0.0003
CADMIUM °	mg/L Cd	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.005*, ≤ 0.05**	0.002
COPPER °	mg/L Cu	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.1	0.002
IRON °	mg/L Fe	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	3.60	-	0.005
LEAD °	mg/L Pb	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.05	0.003
MANGANESE °	mg/L Mn	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	0.569	≤ 1.0	0.002
MERCURY °	mg/L Hg	IN-HOUSE METHOD: UAE.TP.HEM.002 (COLD VAPOUR ATOMIC ABSORPTION SPECTROMETRIC METHOD); SM: 3112 B	ND	≤ 0.002	0.0001
NICKEL °	mg/L Ni	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 0.1	0.005
TOTAL CHROMIUM °	mg/L Cr	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	-	0.005



PARAMETER	UNIT	METHOD OF ANALYSIS	RESULT	REGULATORY STANDARD	DETECTION LIMIT
			BRN-5-SW2 T22AG233-0002		
ZINC ^c	mg/L Zn	IN-HOUSE METHOD: UAE.TP.SW.01 (NITRIC ACID DIGESTION AND DIRECT AIR ACETYLENE FLAME METHOD); SM: 3030 E AND 3111 B	ND	≤ 1.0	0.003
SAMPLE CONDITION WATER'S COLOUR/TURBID SEDIMENT			YELLOW/TURBID YELLOW		

^a : ISO/IEC 17025 ACCREDITED BY THAI INDUSTRIAL STANDARDS INSTITUTE (TISI)

^b : ISO/IEC 17025 ACCREDITED BY DEPARTMENT OF SCIENCE SERVICE (DSS)

^c : VERIFIED BY OWN LABORATORY QUALITY SYSTEM, BUT STILL NOT ACCREDITED

IN-HOUSE : BASED ON STANDARD METHODS FOR THE EXAMINATION OF WATER AND WASTEWATER, APHA, AWWA, WEF, 23rd EDITION, 2017.

SM : STANDARD METHODS FOR THE EXAMINATION OF WATER AND WASTEWATER, APHA, AWWA, WEF, 23rd EDITION, 2017.

REGULATORY STANDARD : SURFACE WATER QUALITY STANDARDS CLASS 3, NOTIFICATION OF THE NATIONAL ENVIRONMENT BOARD, NO.8, B.E. 2537 ISSUED UNDER THE ENHANCEMENT AND CONSERVATION OF NATIONAL ENVIRONMENTAL QUALITY ACT, B.E. 2535, PUBLISHED IN THE ROYAL GOVERNMENT GAZETTE, VOL. 111, PART 16, DATED FEBRUARY 24, B.E. 2537 (1994).

CLASS 3 : MEDIUM CLEAN FRESH SURFACE WATER RESOURCES USED FOR

- (1) CONSUMPTION, BUT PASSING THROUGH ON ORDINARY TREATMENT PROCESS BEFORE USING
- (2) AGRICULTURE

≤ 0.005* : WHEN WATER HARDNESS NOT MORE THAN 100 mg/L AS CaCO₃

≤ 0.05** : WHEN WATER HARDNESS MORE THAN 100 mg/L AS CaCO₃

ND : NON-DETECTABLE

Benjawan V.

(MISS BENJAWAN VIRIYOTHAI)
LABORATORY SUPERVISOR

APRIL 22, 2022

ภาคผนวก ฎ
เอกสารสอบเทียบเครื่องมือ



ใบรับรองสอบเทียบเครื่องมือประจำห้องปฏิบัติการ สำหรับตรวจวัดคุณภาพสิ่งแวดล้อม

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*	Remark
เครื่องมือประจำห้องปฏิบัติการวิเคราะห์ สำหรับวิเคราะห์คุณภาพน้ำใต้ดิน									
1	Inductively Coupled Plasma-Optical Emission Spectrometer (ICP-OES)	กลุ่มโลหะหนัก : ตะกั่ว (Pb), นิกเกิล (Ni),ปรอททั้งหมด (Total Hg), ปรอท (Hg)	Agilent Technologies	System ID:G8015A G8015AA / MY18030001	Agilent Technologies (Thailand) Co.,Ltd.	Preventive Maintenance Checklist	9 Dec 21	8 Dec 22	-
2	Atomic Absorption Spectrometer (AAS)	ทองแดง (Cu),แมงกานีส (Mn),สังกะสี (Zn), เหล็ก (Fe), สารหนู (As), แคดเมียม (Cd),โครเมียมเฮกซะวาเลนท์ (Cr6+)	Agilent Technologies	System ID:G8432A AA240FS / MY13160001	Thailand Institute of Scientific and Technological Research (TISTR).	MTC.ACL. No. 335/64	28 Jan 22	27 Jan 23	-
3	Conductivity Meter	การนำไฟฟ้า(EC)	SI Analytics	Lab955 / 16300356	SPC Calibration Center Co.,Ltd.	C24210091	22 Mar 22	28 Mar 22	-
4	pH Meter	ค่าความเป็นกรด-ด่าง (pH) อุณหภูมิ (Temperature)	Mettler-Toledo	Seven Easy S20 / 1231155210	National Food Institute, Ministry of Industry, Thailand	2103189-002-01	1 Jun 21	13 Jun 22	-
5	pH Meter		Hanna Instrument	HI2211 / 8165345	National Food Institute, Ministry of Industry, Thailand	2102015-001-01	16 Mar 22	16 Mar 22	-
6	Analytical Balance (Repeatability 0.1 mg)	ปิโตรเลียมไฮโดรคาร์บอนทั้งหมด (TPH)	Mettler-Toledo	AB-204S/FACT / 1129361010	National Food Institute, Ministry of Industry, Thailand	2103270-001-01	11 Jun 21	10 Jun 22	-
7	UV-VIS Spectrophotometer	ซัลเฟต (SO ₄ ²⁻)	Agilent Technologies	Cary60 G6860A / MY15410009	DQE Services Co.,Ltd.	SP21-015	29 May 21	28 May 22	-
8	UV-VIS Spectrophotometer		Hitachi	U-1900 / 2021-064	DQE Services Co.,Ltd.	SP21-008	20 Jan 22	19 Jan 23	-
9	UV-VIS Spectrophotometer		Hitachi	U-2900 / 21E22-009	DQE Services Co.,Ltd.	SP22-008	20 Jan 22	19 Jan 23	-

Due Date of Calibration* : กำหนดตามแผนการสอบเทียบประจำปี อย่างน้อยปีละ 1 ครั้ง

Agilent 5110 and 5100 ICP-OES Preventive Maintenance Checklist

Agilent Preventive Maintenance provides factory recommended service for your analytical systems to assure reliable operation and the accuracy of your results. Delivered by highly-trained and certified service engineers using genuine Agilent parts and supplies, Agilent Preventive Maintenance provides everything you need to reduce unplanned downtime and keep your systems operating at their peak.

For more information about Agilent Technologies services please visit our web site using the following URL <http://www.agilent.com/en-us/services/analytical-instrument-services>

Customer Information

- Customers should provide all necessary operating supplies upon request of the engineer.
- For customers using HF applications, the instrument should be returned to its standard sample introduction system.
- A customer representative should be available to the engineer while performing the preventive maintenance procedures.
- Any parts, not included in the Parts Lists section of this document, are not part of the recommended Preventive Maintenance service, nor are they included in the price of this service.
- If a system requires the use of additional or special procedures and/or parts for the instrument service, then these must be ordered separately and charged as a repair, which may incur additional

Service Engineer's Responsibilities

- Only complete/printout pages that relate to the system being serviced.
- Complete empty fields with the relevant information.
- Complete the relevant checkboxes in the checklist using a "X" or tick mark "✓" in the checkbox.
- Complete Not Applicable check boxes to indicate services not delivered, as needed.
- Complete the PM service in the order of the tasks listed.
- Complete the Service Review section together with the customer.

Issued: 4 March 2021, Revision: A.01

Copyright © 2021
Page 1 of 8

Agilent Technologies

เอกสารไม่ควบคุม

Agilent 5110 and 5100 ICP-OES Preventive Maintenance Checklist

System Information

Instrument system name and ID	ICP-OES 5110 VDV
Instrument system site and location	UAE Consultant
List system component product numbers	List the serial numbers of each component
1. 83015A	1. MY 19030001
2.	2.
3.	3.
4.	4.
5.	5.
6.	6.
7.	7.
8.	8.
9.	9.
10.	10.

ICP-OES Configuration table	Circle the type or write in the type if other
Nebulizer Type	SeaSpray OneNeb other
Spray Chamber	Cyclonic Single Pass Cyclonic Double Pass other
Torch	Radial Dual View other
Injector Diameter	2.4mm 1.8mm 1.4mm 0.8mm other
Injector Material	Quartz Ceramic other

Issued: 4 March 2021, Revision: A.01

Copyright © 2021
Page 2 of 8

Agilent Technologies

เอกสารไม่ควบคุม

Agilent 5110 and 5100 ICP-OES Preventive Maintenance Checklist

General Preparation

- Discuss any specific questions or issues with the customer prior to starting.
- Review the instrument logbook.
- Perform general external inspection of system for cleanliness.
- Check for proper installation of safety-related parts, assemblies, sensors etc.
- Check for required firmware/software updates and verify with customers if they would like it installed.
- For HF application systems, if standard sample introduction system was not installed, ask the customer to install it. NA
- Run Instrument Performance test and record results in Instrument Performance Test Results Table - Pre PM.

Inspect and clean the system

- Look for any obvious external damage or problems.
- Inspect water cooling hoses, gas lines and power cord for excessive wear or damage.
- Perform a general internal inspection of the system for excessive dust accumulation, clean if necessary.
- Inspect sample introduction components and record any required maintenance in the Service Engineer Comments and notify the customer as the required actions required.
- Record the instrument operating conditions in the ICP-OES Status Results Table.
- Replace the polychromator purge filter.
- Replace the radial pre-optics window
- Replace the axial pre-optics window for SVDV and VDV instruments.
- Check exhaust flow for the correct positive extraction at the exhaust duct to insure they meet minimum specifications.
- Replace air inlet dust filter.
- Replace high capacity air inlet dust filter element if installed. NA
- Remove and clean instrument water inlet filter.

Agilent Water Recirculator

- Section NOT Applicable
- Drain cooling fluid and remove any particles from the chiller reservoir
- Remove, clean, and reinstall water inlet metal mesh filter if present.
- Re fill with Polyclear Plus cooling fluid.
- Clean the cooling system Air filter and the condenser.

Issued: 4 March 2021, Revision: A.01

Copyright © 2021
Page 3 of 8

Agilent Technologies

เอกสารไม่ควบคุม

Agilent 5110 and 5100 ICP-OES Preventive Maintenance Checklist

SPS 3 Auto Sampler

- Section NOT Applicable
- Power cycle the autosampler and verify successful initialization.
- Inspect X and Z axis belts for wear. Replace is necessary.
- Clean X and Z axis slide shafts.
- Using customer's racks and the Agilent software move the sample probe to the 4 outermost corners and rinse port, ensure that the probe is approximately centered in the vial.

SPS 4 Auto Sampler

- Section NOT Applicable
- Clean the spill tray, rack location mat, end frames and chassis with a damp soft cloth and diluted mild detergent.
- Clean the auto sampler cover panels, if cover kit is installed, with domestic window cleaner
- Check the X-axis and Z-axis drive belts for cracks, splits, damaged teeth, excessive fraying, color changes or degradation from fumes.
- Check the X-axis, Theta-axis and Z-axis FFC cables for cracks, incorrect positioning, damaged edges or damaged connectors.
- Pump Tubing Replacement. Replace peristaltic pump tubing. Replace all tubing that goes from the rinse station to the pump and from the pump to the waste/rinse bottles

AVS 4, 6, 7

- Section NOT Applicable
- Replace valve rotor seal
- Check fittings for signs of leaks
- Check tubing including autosampler tubing for kinks or excessive wear
- Check high flow pump for signs of leaks

Instrument Adjustment

- Check position of Zn peak, adjust if required.
- Check Argon Ratio, adjust to specified value if required.
- Perform Detector Calibration.
- Perform Instrument Calibration.
- Run Instrument Performance Test and record results in Instrument Performance Test Results Table - Post PM.
- For systems using ICP Expert version 7.3 and above run the following Instrument tests and record the result in the Instrument Test Results Table
 - Subsystem Communications Test
 - Air Flow

Issued: 4 March 2021, Revision: A.01

Copyright © 2021
Page 4 of 8

Agilent Technologies

เอกสารไม่ควบคุม

**Agilent 5110 and 5100 ICP-OES
Preventive Maintenance Checklist**

- ☒ Water Flow
- ☒ Gas Flows
- ☒ RF Generator
- ☒ Camera Test
- ☒ Optics Test
- ☒ Nebulizer Test

Instrument Performance Test Results Table

Note: These measurements do not form part of any specification and are for reference only.

	Pre PM Sensitivity Check		Post PM Sensitivity Check	
	Radial	Axial *	Radial	Axial*
Zn 213.857 nm SRBR	4092.3	7954.1	4192.8	7920.7
Mn 257.610 nm SRBR	11415.2	30894.7	11493.6	34460.9
Al 396.152 nm SBR	7.8	15.7	8.7	13.5
K 766.491 nm SBR	5.3	39.9	5.7	44.6

* Axial result is not applicable for G8016AA, G8012AA Radial View instruments.

Instrument Test Results Table

Note: The Instrument Test results are for systems using ICP Expert version 7.3 and above only.

Instrument Test	Result
Subsystem Communications Test	PASS
Air Flow	PASS
Water Flow	PASS
Gas Flows	PASS
RF Generator	PASS
Camera Test	PASS
Optics Test	PASS
Nebulizer test	PASS

Issued: 4 March 2021, Revision: A.01

Copyright © 2021
Page 5 of 8

Agilent Technologies

เอกสารไม่ควบคุม

**Agilent 5110 and 5100 ICP-OES
Preventive Maintenance Checklist**
ICP-OES Status Results Table

Note: These measurements do not form part of any specification and are for reference only.

Measurement	Standby Mode	Plasma On
Mains Voltage	115.717 VAC	194.510 VAC
Mains Current	0.219 A	3.233 A
Instrument Temperature	23.4 °C	23.5 °C
RF Air Flow (sensor speed)	14.0 Hz	19.0 Hz
Plasma Exhaust Temperature	No measurement	65.0 °C
Water Flow Oscillator	No measurement	2.03 L/min
Water Flow Detector	0.00 L/min	1.37 L/min
Water Inlet Temperature	19.2 °C	12.8 °C
Polychromator Temperature	35.0 °C	35.0 °C
CCD Temperature	26.9 °C	-39.7 °C
Thermal Stabilizer	35.0 °C	35.0 °C
Argon Supply Pressure	614.25 kPa	629.97 kPa
Purge Gas Supply Pressure*1	641.74 kPa	654.67 kPa
Option Gas Supply Pressure*1	- kPa	- kPa
Nebulizer Flow	No measurement	0.90 L/min
Nebulizer Back Pressure	No measurement	163.05 kPa
Plasma Gas Flow	No measurement	15.00 L/min
Auxiliary Gas Flow	No measurement	1.20 L/min
RF Power	No measurement	1207.1 W
RF Supply Current	No measurement	3.233 A
RF Supply Voltage	No measurement	194.510 V

*1 If option installed

Issued: 4 March 2021, Revision: A.01

Copyright © 2021
Page 6 of 8

Agilent Technologies

เอกสารไม่ควบคุม

**Agilent 5110 and 5100 ICP-OES
Preventive Maintenance Checklist**
ICP-OES Parts List Table

Part description	Part Number	Product / Model # where used	Quantity Consumed
Axial Pre-Optic Window	G8010-68014	G8010A, G8011A, G8014A/G8015A	1
Radial Pre-Optic Window	G8010-68015	All	1
Polyclear Plus Cooling Fluid	G3292-80012	Agilent Water Recirculator	-
Purge Gas Filter	G8010-60136	All	1
Air inlet filter	G8000-68002	All	1
High Capacity Air Filter	G8010-60189	Optional	-
Rotor seal for 6-7 port valve for AVS6/7	G8494-60002	G8494A/G8495	-
Rotor seal for 4 port valve for AVS4	G8493-60002	G8493A	-
Rinse solution to rinse station 2.5mm id x 1m	G8410-80123	SPS 4	-
Barb connector 2.5mm-1.5mm ID	G8410-80124	SPS 4	-
PVC waste tubing, 8mm od x 5mm id, 2m	G8410-80122	SPS 4	-
Additional Parts may be required from engineers stock:			
X axis drive belt	5410047500	SPS 3	-
Z axis drive belt	5410047400	SPS 3	-
Peristaltic pump tubing, PVC SolvaFlex, 3 bridged,	3710049000	SPS 4	-

Restore system

For HF applications, ask the customer to reinstall their sample introduction system.

Leave system in an idle state: on and purging.

Guidance: If the PM service is performed prior to a qualification service, then use the qualification procedure as a guide for final instrument set up and checkout.

Service Review

- ☒ Affix the PM sticker to the system or instrument logbook based on the customer's request.
- ☒ Complete the Service Engineer Comments section below if there are additional comments.

Issued: 4 March 2021, Revision: A.01

Copyright © 2021
Page 7 of 8

Agilent Technologies

เอกสารไม่ควบคุม

**Agilent 5110 and 5100 ICP-OES
Preventive Maintenance Checklist**

- ☒ Review the service and any test results with the customer.
- ☒ If the Instrument firmware was updated, record the details of the change in the Service Engineer's Comments box below or if necessary, in the customer's IQ records.

Service Engineer Comments (optional)

If there are any specific points you wish to note as part of performing the installation or other items of interest for the customer, please write in this box.

Other Important Customer Web Links

How to get information on your product:

- ☐ Literature Library - <http://www.agilent.com/en-us/products/icp-oes/icp-oes-systems/5110-icp-oes#literature>
- ☐ Need to know more? - <http://www.agilent.com/crosslab/university/>
- ☐ Need technical support, FAQs? - <http://www.agilent.com/en-us/support/landing/icp-oes>
- ☐ Need supplies? - www.agilent.com/chem/supplies

Service Completion

Service request number 600439717 Date service completed 09/12/21

Agilent signature Nelson L. Customer signature Aphorn Onkong

Document part number: G8014-90075

Issued: 4 March 2021, Revision: A.01

Copyright © 2021
Page 8 of 8

Agilent Technologies

เอกสารไม่ควบคุม

Report Summary	
Instrument Model	Agilent 5100/5110 VDV ICP-OES
Instrument ID	G8011A/G8015A
Instrument Serial Number	MY18030001
Software Version	7.3.1.9507
Firmware Version	3442
Tested By	Nukoon L.
Test Completed On	12/9/2021 9:14:59 AM
Result Summary	
Subsystem Communications Test	Skipped
Air Flow Test	Skipped
Water Flow Test	Skipped
Gas Flows Test	Skipped
RF Generator Test	Skipped
Camera Test	Skipped
Optics Test	Skipped
Advanced Valve System Test	Skipped
Resolution Test	Pass
Sensitivity Test	Pass
Precision Test	Pass

เอกสารไม่ควบคุม

Resolution Test			Pass
Element Wavelength	Specification	Width	
N (174.213 nm)	≤ 9.40	7.27	
As (188.980 nm)	≤ 8.20	8.23	
C (193.027 nm)	≤ 11.50	8.26	
Mo (202.032 nm)	≤ 8.20	6.42	
Cr (206.158 nm)	≤ 13.40	9.27	
Zn (213.857 nm)	≤ 8.70	6.77	
Pb (220.353 nm)	≤ 9.50	7.12	
Co (228.615 nm)	≤ 17.20	11.88	
Ba (230.424 nm)	≤ 9.40	7.36	
Mn (257.610 nm)	≤ 13.30	9.52	
Mn (260.568 nm)	≤ 20.30	14.30	
Cr (267.716 nm)	≤ 11.00	7.99	
Cu (324.754 nm)	≤ 25.00	19.06	
Cu (327.395 nm)	≤ 14.20	11.32	
Sr (338.071 nm)	≤ 33.50	24.39	
Ba (455.403 nm)	≤ 44.00	33.86	
Sr (460.733 nm)	≤ 36.00	17.38	
Ba (493.408 nm)	≤ 36.00	25.53	
Ba (614.171 nm)	≤ 42.00	24.99	
Ar (675.283 nm)	≤ 74.00	59.49	
K (766.491 nm)	≤ 80.00	65.27	

เอกสารไม่ควบคุม

Sensitivity Test						Pass
Radial						
Element Wavelength	Specification	Method	Ratio	Standard	Blank	
As (188.980 nm)	≥ 46.0	SRBR	167.2	1131.3	42.4	
Se (196.026 nm)	≥ 41.0	SRBR	119.1	1177.1	84.2	
Zn (213.857 nm)	≥ 1421.0	SRBR	4082.3	49908.2	148.6	
Pb (220.353 nm)	≥ 46.0	SRBR	191.1	2682.8	172.6	
Mn (257.610 nm)	≥ 3518.0	SRBR	11415.2	265002.2	536.8	
Al (396.152 nm)	≥ 3.4	SBR	7.8	49838.0	5676.5	
Ba (493.408 nm)	≥ 34.0	SBR	116.1	1999041.4	17066.5	
K (766.491 nm)	≥ 1.8	SBR	5.3	101078.4	16104.6	
Axial						
Element Wavelength	Specification	Method	Ratio	Standard	Blank	
As (188.980 nm)	≥ 208.0	SRBR	252.9	3214.2	147.0	
Se (196.026 nm)	≥ 159.0	SRBR	216.2	3839.7	272.2	
Zn (206.200 nm)	≥ 234.0	SRBR	1203.3	14046.1	133.7	
Zn (213.857 nm)	≥ 1743.0	SRBR	7856.1	171323.1	472.9	
Cd (214.439 nm)	≥ 4227.0	SRBR	7054.9	129539.3	335.4	
Pb (220.353 nm)	≥ 320.0	SRBR	531.7	13218.2	566.2	
Mn (257.610 nm)	≥ 10625.0	SRBR	30884.7	1314844.0	1807.4	
Cr (267.716 nm)	≥ 1048.0	SRBR	4442.1	174420.3	1515.1	
Cu (324.754 nm)	≥ 19.0	SBR	50.7	374603.6	7249.0	
Al (396.152 nm)	≥ 6.0	SBR	15.7	279915.3	16790.4	
Ba (493.408 nm)	≥ 60.0	SBR	209.7	10899956.6	51728.3	
K (766.491 nm)	≥ 24.0	SBR	38.9	1983197.5	49746.6	

เอกสารไม่ควบคุม

Precision Test			Pass
Radial			
Element Wavelength	Specification	Measured Value % RSD	
As (188.980 nm)	≤ 2.60	0.81	
Se (196.026 nm)	≤ 2.60	1.21	
Zn (213.857 nm)	≤ 1.50	0.39	
Pb (220.353 nm)	≤ 2.60	0.41	
Mn (257.610 nm)	≤ 1.50	0.45	
Al (396.152 nm)	≤ 1.50	0.41	
Ba (493.408 nm)	≤ 1.50	0.51	
K (766.491 nm)	≤ 1.50	0.36	
Axial			
Element Wavelength	Specification	Measured Value % RSD	
As (188.980 nm)	≤ 1.50	0.51	
Se (196.026 nm)	≤ 1.50	0.73	
Zn (206.200 nm)	≤ 1.50	0.30	
Zn (213.857 nm)	≤ 1.50	0.37	
Cd (214.439 nm)	≤ 1.50	0.36	
Pb (220.353 nm)	≤ 1.50	0.28	
Mn (257.610 nm)	≤ 1.50	0.63	
Cr (267.716 nm)	≤ 1.50	0.30	
Cu (324.754 nm)	≤ 1.50	0.54	
Al (396.152 nm)	≤ 1.50	0.45	
Ba (493.408 nm)	≤ 1.50	0.64	
K (766.491 nm)	≤ 1.50	0.56	

เอกสารไม่ควบคุม

Report Summary

Instrument Model	Agilent 5100/5110 VDV ICP-OES
Instrument ID	G8011A/G8015A
Instrument Serial Number	MY18030001
Software Version	7.3.1.9507
Firmware Version	3442
Tested By	Nukoon L.
Test Completed On	12/9/2021 12:55:49 PM

Result Summary

Subsystem Communications Test	Skipped
Air Flow Test	Skipped
Water Flow Test	Skipped
Gas Flows Test	Skipped
RF Generator Test	Skipped
Camera Test	Skipped
Optics Test	Pass
Advanced Valve System Test	Skipped
Resolution Test	Pass
Sensitivity Test	Pass
Precision Test	Pass

Optics Test

Pass

	Radial	Axial
Intensity	5296135	5755042
Wavelength	737.212	737.212

Page 1 of 4

เอกสารไม่ควบคุม

Resolution Test

Pass

Element Wavelength	Specification	Width
N (174.213 nm)	≤ 9.40	7.22
As (188.980 nm)	≤ 8.20	6.15
C (193.027 nm)	≤ 11.50	8.22
Mo (202.032 nm)	≤ 8.20	6.33
Cr (206.158 nm)	≤ 13.40	9.21
Zn (213.857 nm)	≤ 8.70	6.87
Pb (220.353 nm)	≤ 9.50	7.02
Co (228.615 nm)	≤ 17.20	11.81
Ba (230.424 nm)	≤ 9.40	7.46
Mn (257.610 nm)	≤ 13.30	9.49
Mn (260.568 nm)	≤ 20.30	14.19
Cr (267.716 nm)	≤ 11.00	7.90
Cu (324.754 nm)	≤ 25.00	18.92
Cu (327.395 nm)	≤ 14.20	11.32
Sr (338.071 nm)	≤ 33.50	24.29
Ba (455.403 nm)	≤ 44.00	33.68
Sr (460.733 nm)	≤ 36.00	17.64
Ba (493.408 nm)	≤ 36.00	25.56
Ba (614.171 nm)	≤ 42.00	24.75
Ar (675.283 nm)	≤ 74.00	59.18
K (766.491 nm)	≤ 80.00	65.19

Page 2 of 4

เอกสารไม่ควบคุม

Sensitivity Test

Pass

Radial	Element Wavelength	Specification	Method	Ratio	Standard	Blank
	As (188.980 nm)	≥ 46.0	SRBR	154.8	1242.3	58.5
	Se (196.026 nm)	≥ 41.0	SRBR	117.4	1259.6	97.9
	Zn (213.857 nm)	≥ 1421.0	SRBR	4192.8	52402.6	155.3
	Pb (220.353 nm)	≥ 46.0	SRBR	196.4	2814.2	179.9
	Mn (257.610 nm)	≥ 3518.0	SRBR	11993.6	281210.1	547.6
	Al (396.152 nm)	≥ 3.4	SBR	8.7	55103.6	5662.9
	Ba (493.408 nm)	≥ 34.0	SBR	125.4	2152916.9	17032.2
	K (766.491 nm)	≥ 1.8	SBR	5.7	107906.7	16079.8

Axial

Element Wavelength	Specification	Method	Ratio	Standard	Blank
As (188.980 nm)	≥ 208.0	SRBR	297.5	4054.8	170.4
Se (196.026 nm)	≥ 159.0	SRBR	260.2	4794.9	298.5
Zn (206.200 nm)	≥ 234.0	SRBR	1305.9	16162.3	150.3
Zn (213.857 nm)	≥ 1743.0	SRBR	8920.7	200915.6	504.7
Cd (214.439 nm)	≥ 4227.0	SRBR	7958.3	149327.5	350.4
Pb (220.353 nm)	≥ 320.0	SRBR	606.7	15244.5	584.0
Mn (257.610 nm)	≥ 10625.0	SRBR	34460.9	1493092.8	1872.5
Cr (267.716 nm)	≥ 1048.0	SRBR	5018.6	198000.6	1532.6
Cu (324.754 nm)	≥ 19.0	SBR	57.5	423683.7	7248.6
Al (396.152 nm)	≥ 6.0	SBR	18.5	320004.9	16441.4
Ba (493.408 nm)	≥ 60.0	SBR	233.3	11882915.4	50714.5
K (766.491 nm)	≥ 24.0	SBR	44.6	2218974.4	48657.9

Page 3 of 4

เอกสารไม่ควบคุม

Precision Test

Pass

Radial	Element Wavelength	Specification	Measured Value % RSD
	As (188.980 nm)	≤ 2.60	1.38
	Se (196.026 nm)	≤ 2.60	0.91
	Zn (213.857 nm)	≤ 1.50	0.38
	Pb (220.353 nm)	≤ 2.60	0.44
	Mn (257.610 nm)	≤ 1.50	0.43
	Al (396.152 nm)	≤ 1.50	0.38
	Ba (493.408 nm)	≤ 1.50	0.66
	K (766.491 nm)	≤ 1.50	0.36

Axial

Element Wavelength	Specification	Measured Value % RSD
As (188.980 nm)	≤ 1.50	0.61
Se (196.026 nm)	≤ 1.50	0.52
Zn (206.200 nm)	≤ 1.50	0.36
Zn (213.857 nm)	≤ 1.50	0.33
Cd (214.439 nm)	≤ 1.50	0.41
Pb (220.353 nm)	≤ 1.50	0.36
Mn (257.610 nm)	≤ 1.50	0.74
Cr (267.716 nm)	≤ 1.50	0.25
Cu (324.754 nm)	≤ 1.50	0.71
Al (396.152 nm)	≤ 1.50	0.44
Ba (493.408 nm)	≤ 1.50	0.73
K (766.491 nm)	≤ 1.50	0.97

Page 4 of 4

เอกสารไม่ควบคุม

Report Summary		
Instrument Model	Agilent 5100/5110 VDV ICP-OES	
Instrument ID	G8011A/G8015A	
Instrument Serial Number	MY18030001	
Software Version	7.3.1.9507	
Firmware Version	3442	
Tested By	Nukoon L.	
Test Completed On	12/9/2021 1:34:10 PM	
Result Summary		
Subsystem Communications Test	Pass	
Air Flow Test	Pass	
Water Flow Test	Pass	
Gas Flows Test	Pass	
RF Generator Test	Pass	
Camera Test	Pass	
Optics Test	Skipped	
Advanced Valve System Test	Skipped	
Resolution Test	Skipped	
Sensitivity Test	Skipped	
Precision Test	Skipped	
Subsystem Communications Test	Pass	
Air Flow Test	Pass	
30% Air Flow (relative speed)	75% Air Flow (relative speed)	
15.00	19.00	
Water Flow Test	Pass	
RF Water Flow(L/min)	Camera Water Flow (L/min)	Water Inlet Temperature (°C)
1.98	1.36	17.16

Page 1 of 2

เอกสารไม่ควบคุม

Gas Flows Test			Pass		
Nebulizer Target Flow	Actual Flow	Back Pressure	Auxiliary Target Flow	Actual Flow	Back Pressure
0.70	0.70	203.80	2.00	1.99	108.66
Makeup Target Flow	Actual Flow	Back Pressure	Plasma Target Flow	Actual Flow	Back Pressure
2.00	2.00	113.89	18.00	17.93	24.24
RF Generator Test			Pass		
RF Power Supply Test	Passed				
RF Power Supply (V)	141.475				
RF Oscillator Test	Passed				
RF Oscillator Frequency (MHz)	25.874				
Work Coil Current (A)	45.931				
RF Power Supply Current (A)	2.000				
Camera Test			Pass		
	Integration Time (ms)	Standard Deviation	Status		
Electronic Offset Test	1000	5.261	Passed		
Dark Current Test	6000	0.734	Passed		
Array Test	5	0.024	Passed		
Linearity Test		0.118	Passed		

Page 2 of 2

เอกสารไม่ควบคุม



Agilent 55 240 280 Series Atomic Absorption Spectroscopy Systems

Preventive Maintenance Checklist

Agilent Preventive Maintenance provides factory recommended service for your analytical systems to assure reliable operation and the accuracy of your results.

Delivered by highly trained and certified service engineers using genuine Agilent parts and supplies, Agilent Preventive Maintenance provides everything you need to reduce unplanned downtime and keep your systems operating at their peak. This checklist will be completed at the end of the service and provided to you as a record of the installation.

Note: While non-current production AA instrument and/or accessory models are not covered specifically in this document it can be used as a basic reference.

For more information about Agilent Technologies services please visit our web site using the following URL <http://www.agilent.com/en-us/services>

Introduction

Customer Information

- Customers should provide all necessary operating supplies upon request of the engineer.
- A customer representative should be available to the engineer while performing the preventive maintenance procedures.
- Any parts, not included in the Parts Lists section of this document, are not part of the recommended Preventive Maintenance service, nor are they included in the price of this service.
- If a system requires the use of extra or special procedures and/or parts for the maintenance service, then these must be ordered separately and charged as a repair, which may incur additional costs.



Instrument Preventive Maintenance Checklist

Important Customer Web Links

- For more information about Agilent Technologies services, please visit our website using the following URL: <http://www.agilent.com/en-us/products/crosslab-instrument-services/service-repair>
- To access Agilent University, visit <http://www.agilent.com/crosslab/university/> to learn about training options, which include online, classroom and onsite delivery. A training specialist can work directly with you to help determine your best options.
- A useful Agilent Resource Center web page is available, which includes short videos on maintenance, quick lists of consumables for new instruments, and other valuable information. Check out the Resource Page here: <https://www.agilent.com/en-us/agilentresources>
- Need technical support, FAQs, supplies? – visit our Support Home page at <http://www.agilent.com/search/support>
- Get answers. Share insights. Build connections: Join the Agilent Community at <https://community.agilent.com/welcome>

Service Engineer's Responsibilities

- Contact the customer and ensure that all necessary supplies are available before the preventive maintenance visit.
 - Confirm the ability of the instrument to deliver continued safe operation as established via the Agilent AA safe operation flow chart. (Refer directly to the AA 55/240/280 Preventive Maintenance Scope of Work to make this decision.)
 - Only select those pages that relate to the system or module being serviced.
 - Complete empty fields with the relevant information.
 - Complete the relevant checkboxes in the checklist using either a "X" or tick mark "✓".
 - Check "Section not applicable" check boxes to indicate services/tasks not delivered, as appropriate.
 - Complete the Preventive Maintenance service in the order of the tasks listed.
 - Complete the Service Review section together with the customer.
 - Complete the fields for page numbers at the foot of each selected page
 - Complete the total number of pages field in the Service Completion section
 - Ask the customer to sign the Service Completion section including the customer's and your signature.
- This information is subject to change without notice.



เอกสารไม่ควบคุม



เอกสารไม่ควบคุม

Instrument Maintenance

System Information

- ☐ Check this box if an instrument configuration report is attached instead of completing the table.

Instrument System Name and ID	240FS AAS
Instrument System Site and Location	United Analyst and Engineering Consultant

List System Component Product Numbers	List the Serial Numbers of each Component
1. G 2432 A	MY 13160001
2.	
3.	
4.	
5.	
6.	
7.	
8.	
9.	

Preparation, Safe operation and Initial performance checks

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021



- ☐ Agilent AA safe operation flow chart inspections (to determine if the PM can be performed).

NOTE: If by following the flow chart the instrument is deemed to be unsafe for continued use you MUST NOT continue PM work. Inform the customer immediately of the Agilent recommendation that use of the instrument be discontinued.

- ☒ Discuss any specific issues with the customer before starting.
- ☐ For HF application systems, if standard sample introduction system was not installed, ask the customer to install it. **N/A**
- ☒ Review the instrument logbook for recorded problems and comments.
- ☒ Save instrument control settings before starting the procedure.
- ☒ Perform a general inspection of the system for cleanliness.
- ☒ Check for proper installation of parts, assemblies, sensors etc.
- ☒ Check system for required installation of components, settings as defined by current Service Notes
- ☒ Check for required firmware updates and verify with customers if they would like them installed.
- ☒ Use SVD to perform a Full Wavelength Scan for Cu HCL - "As found test_1"
- ☒ Perform a Basic Cu ABS test - "As found test_2"
- ☒ Print the Details page or screen captures of the test results and attach to the end of this checklist.

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021



Preventive Maintenance Procedures

FLAME SYSTEM section

- ☐ Section not applicable

Electronic components

- ☒ Review and confirm instrument configuration data in SVD
- ☒ Confirm power supply voltages using the **SVD Power Supply diagnostic**.
- ☒ For Dual Beam instruments - Confirm RBC frequency using the **SVD RBC frequency diagnostic**.

Mechanical components

- ☒ Check the burner adjuster controls for complete and free movement. If the burner adjuster needs lubrication, use Molykote 321 or mineral-based molybdenum disulphide grease.
- ☒ Run SVD tests to exercise all motor drives over the full range of their travel:
 - ☒ Monochromator drive
 - ☒ Slit drive
 - ☒ Lamp selector
 - ☐ ABA **N/A**

Optics components

- ☒ Check that external optical surfaces are clean - Clean or replace as required.
- ☒ Use SVD and perform **Mono Wavelength Correction**.
- ☒ Use SVD and perform **Slit Calibration**.
- ☒ Use SVD and perform **Grating Squareness Diagnostic**.
- ☒ Use SVD and perform **Zero Order Offset/Mono Correction**.
- ☒ Use SVD and perform **Wavelength Repeatability**.
- ☒ Physically inspect selected HC lamps (customer to supply per their choice) and measure the % Gain for each lamp. Advise customer if lamps are showing emission degradation due to age.
- ☒ Check that the signal energy of the D2 and HC lamps track properly. Advise customer if their D2 lamp is showing emission degradation due to age.

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021



Sample Introduction and Atomization

- ☒ Inspect the burner interlock plate to ensure that the interlock pin is secure and correct for the burner type.
- ☒ Clean the burner slot with a clean white card.
- ☒ Check the uniformity of the slot width.
- ☒ Clean the burner if required.
- ☒ Change the burner o-ring.
- ☒ Clean the nebulizer, spray chamber and liquid trap.
- ☒ Change all o-rings and seals in the nebulizer, nebulizer block and spray chamber.
- ☒ Check that the pressure relief bung releases readily.
- ☒ Change o-rings on the fuel and oxidant delivery bars
- ☒ Leave the liquid trap EMPTY and verify the flame will not ignite in this state.
- ☒ Refill liquid trap and check that overflow drains freely into the drain/waste tube.
- ☒ Check the drain/waste tube for good drainage. It should not have tight bends, kinks or loops and the lower end must be above the liquid level in the waste vessel
- ☒ Check and clean the igniter electrode

Gas handling components and safety interlocks

- ☒ Pressure test for leaks
- ☒ Leak test gasbox internal components and connections
- ☒ Check safety interlock status and operation using the **SVD interlock monitoring diagnostic**.

Analytical performance for Flame systems

- ☒ Ignite a flame.
- ☒ Check that you can adjust the nebulizer uptake rate from 4 to 6.5 mL per minute.
- ☒ Optimize the instrument ready to perform Cu sensitivity test.
- ☒ Create a manual method to perform a Basic Cu ABS test - "Final Performance Testing"
- ☒ Run a PM completed sensitivity test for a 5 ppm copper sample and record the results in the AA PM Performance test results and measurements table.

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021



FURNACE SYSTEM section

☒ Section not applicable**Electronic components**

- ☐ Review and confirm instrument configuration data in SVD
- ☐ Confirm power supply voltages using the **SVD Power Supply diagnostic**.

Mechanical components

- ☐ Run SVD tests to exercise all motor drives over the full range of their travel:
 - ☐ Monochromator drive
 - ☐ Slit drive
 - ☐ Lamp selector

Optics components

- ☐ Check that external optical surfaces are clean – Clean or replace as required.
- ☐ Use SVD and perform **Mono Wavelength Correction**.
- ☐ Use SVD and perform **Slit Calibration**.
- ☐ Use SVD and perform **Grating Squareness Diagnostic**.
- ☐ Use SVD and perform **Zero Order Offset/Mono Correction**.
- ☐ Use SVD and perform **Wavelength Repeatability**.
- ☐ Physically inspect selected HC lamps (customer to supply per their choice) and measure the % Gain for each lamp. Advise customer if lamps are showing emission degradation due to age.

Gas handling, water system and workhead component checks

- ☐ Inspect the GTA workhead gas hoses and connections for leaks.
- ☐ Pressure test for gas leaks
- ☐ If the cooler system is accessible (stand-alone) check for correct operation and coolant/water level – this includes any temperature and pressure settings plus filter cleaning (air flow and water).
- ☐ Inspect the GTA workhead water hoses and connections for leaks.
- ☐ Check all graphite components and replace if necessary.

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021



- ☐ Tube
- ☐ Electrodes
- ☐ Shroud

- ☐ Check and clean the end windows on the workhead.
- ☐ Check safety interlock operation.

Analytical performance for Furnace systems

- ☐ Optimize the instrument ready to perform Cu sensitivity test.
- ☐ Run the sensitivity test for a 25 ppb copper sample and record the results in the results table.

PSD autosampler accessory for Furnace systems

- ☐ Section NOT Applicable
- ☐ Check condition of the PSD capillary – replace if necessary
- ☐ Check condition and operation of PSD syringe – ensure it does not have air locks and bubbles.
- ☐ Change PSD rinse bottle o-ring.
- ☐ Check and clean the rinse vessel.
- ☐ Check the drain tube for good drainage. It should not have tight bends, kinks or loops and the lower end must be above the liquid level in the waste vessel.
- ☐ Ensure that the waste vessel is suitable for use with the furnace system.

Sample introduction pump system (SIPS) accessory

- ☐ Section NOT Applicable
- ☐ Re-torque screws securing the hubs, presser arms and pump rotors.
- ☐ Adjust each roller so that it rotates freely.
- ☐ Wipe clean the pump rotor rollers and pump bands with a dry clean cloth.
- ☐ Ensure that the presser arms and the surfaces near the pump are free from dirt and spills.
- ☐ Remove the pump module rear cover and check for the incursion of liquids and any signs of corrosion.
- ☐ Re-torque the nuts that fasten the motor mounting plates to the chassis.
- ☐ Check clips securing the diluents holder and replace if necessary.
- ☐ Disconnect, clean T-piece, and reassemble the tubing using the following steps.

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021



- ☐ Remove the T-piece by disconnecting the pump tubes, the pump bands and all other tubing.
- ☐ Place the T-piece in an ultrasonic bath containing strong detergent 1-5% Decon 30 or similar, for approximately 5-10 minutes.
- ☐ Wash the T-piece under a tap with a strong flow of water.
- ☐ Rinse with distilled water through all of the inlets in the reverse direction to normal sample flow.
- ☐ Reassemble.

Sample preparation system (SPS 4) accessory☒ Section NOT Applicable

The Agilent SPS 4 autosampler is designed to need minimal maintenance.

The following maintenance requirements are suggested to maintain the performance of the autosampler.

- ☐ Cleaning the spill tray, rack location mat, end frames and chassis accessories with a damp soft cloth and diluted mild detergent.
- ☐ Cleaning the autosampler cover panels with domestic window cleaner.
- ☐ Checking the X- axis and Z- axis drive belts for cracks, splits, damaged teeth, excessive fraying, color changes or degradation from fumes.
- ☐ Check the X- axis, Theta- axis and Z- axis FFC cables for cracks, incorrect positioning, damaged edge or damaged connectors.

NOTE: The autosampler requires no extra lubrication throughout its lifetime.

For further details refer to the SPS 4 service manual G8410-90050.

Sample preparation system (SPS 3) accessory☒ Section NOT Applicable

- ☐ Check the x-axis and z-axis timing belts – Replace if there is are any cracks, splits or color deterioration and belt tension.
- ☐ Check belt tensions - adjust if required
- ☐ Check the lubrication pad for single x-axis shaft. If pad is dry or customer has observed any vibration or erratic movements of the x-axis carriage, add 1 mL of Dow Corning 200 8 Fluid, 200 CS into the well.
- ☐ Check the auto-sampler ability to find tube positions - Calibrate if required.
- ☐ Clean the exterior surfaces of the accessory with soft lint free cloth. This cloth can be dampened with warm water or a mild detergent. Do not use organic solvents or abrasive cleaning agents.

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021

**Vapor generation accessory VGA (hydride generator)**

- ☐ Section NOT Applicable
- ☒ Inspect VGA gas supply hose.
- ☒ Inspect/replace VGA pump tubing.
- ☒ Check low gas pressure interlock setting– adjust if required.
- ☒ Check precision orifice gas flow setting – adjust if required.
- ☒ Check gas regulator pressure to 46 psi (325 kPa) – adjust if required.
- ☐ Clean the exterior surfaces of the accessory with soft lint free cloth. This cloth can be dampened with warm water or a mild detergent. Do not use organic solvents or abrasive cleaning agents. **N/A**

UltraAA lamp accessory (external)

- ☒ Section NOT Applicable
- ☐ Check the condition of the power cable.
- ☐ Clean the exterior surfaces of the accessory with soft lint free cloth. This cloth can be dampened with warm water or a mild detergent. Do not use organic solvents or abrasive cleaning agents.

Restore System


- ☐ If you have altered the customer's instrumentation during the course of PM, restore to the original status to allow the customer to conduct their normal activities (e.g., reload the customer's method.)

Guidance

If the PM service is performed prior to a qualification service, then use the qualification procedure as a guide for final instrument set up and checkout.

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021

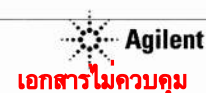

Signature Page

Service Review

- ☒ Attach available reports/printouts of all tests to this documentation.
- ☒ Record the Preventive Maintenance service activity in the customer's records/logbook.
- ☒ Update/reset instrument maintenance counters as appropriate.
- ☒ Affix the PM sticker to the system or instrument logbook based on the customer's request.
- ☒ Complete the Service Engineer Comments section if there are additional comments.
- ☒ Review this service, parts replaced, and test results obtained with the customer.
- ☒ If the instrument firmware was updated, record the details of the change in the Service Engineer's Comments box or if necessary, in the customer's IQ records.

Test Results

Test Description	Expected Test Result	Actual Test Result
Flame optics PMT Gain test		
For copper at 324.8 nm, 4 mA, 0.5 nm slit width	< 55 %	42 %
Flame performance test with 5 ppm copper sample		
Air /acetylene, mixing paddle removed	Abs value > 0.5	0.5963
Air /acetylene, mixing paddle installed, 10 replicates	%RSD < 1.0	0.4
Deuterium furnace optics PMT Gain test		
For copper at 324.8 nm, 4 mA, 0.5 nm slit width	< 55 %	
Deuterium furnace performance test with 25 ppb copper sample (327.4 nm)		
Precision %RSD	≤ 4.0 %	-
Abs value	≥ 0.15	-
Zeeman furnace analytical performance: 25 ppb copper sample (327.4 nm)		
Precision %RSD	≤ 4.0 %	-
Abs value	≥ 0.10	-
MSR%	≥ 70 %	-
As V&A-77 test Abs.	0.2 Abs	0.2156 Abs

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021



AA consumable and parts list table

Part Description	Part Number	Product/Model # where used	PM supplied or Consumable	Instrument-Type
Test Solution - Cu 5ppm solution	6610030100	50 55 140 240 280	PM supplied	Common
Test Solution - Blank solution	5190-7001	50 55 140 240 280	PM supplied	Common
Copper, 1000 ug/ml, 100ml	5190-8279	50 55 140 240 280	*	Common
Kit, Mkr 7 O-rings, aqueous, complete set	9910093400	50 55 140 240 280	PM supplied	Flame
Organic Kit	9910093500	50 55 140 240 280	PM supplied	Flame
Wire Nebulizer Cleaning	9910024700	50 55 140 240 280	consumable	Flame
Tubing-Capillary Std Nebs	9910024800	50 55 140 240 280	consumable	Flame
Capillary Tube Hivac Neb (3) (organics only)	9910044000	50 55 140 240 280	consumable	Flame
Glass impact beads (5/pk)	9910025700	50 55 140 240 280	consumable	Flame
Teflon impact beads (5/pk) (organics only)	9910053300	50 55 140 240 280	consumable	Flame
Burner cleaning strip (100/pk)	9910053900	50 55 140 240 280	consumable	Flame
Window UV silica - round (right side)	2010082600	50 55 140 240 280	PM supplied	Common
Window UV silica - rectangular (left side)	2010082500	50 55 140 240 280	PM supplied	Common
Pad adhesive window (round)	4910012700	50 55 140 240 280	PM supplied	Common
Pad adhesive window (rectangular)	4910012800	50 55 140 240 280	PM supplied	Common
Electrode kit (1 pr) (D2)	6310003400	GTA120	PM supplied	Furnace
Shroud (D2)	6310003100	GTA120	PM supplied	Furnace
Zeeman electrode kit (1 pr)	6310003500	GTA120	PM supplied	Furnace
Zeeman shroud	6310003600	GTA120	PM supplied	Furnace
O-ring PSD rinse bottle	6910025900	PSD120	PM supplied	Furnace

* For engineers who only service AA instruments 5190-8279 can be used as a cheaper alternative for 6610030100.

Items classified as PM supplied in the above table are included in the standard PM Those classified as consumable should be provided by the customer or charged to the customer if supplied by the Agilent service engineer.

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021



SVD Results Report



Report ID: Diagnostic Start Time: 28/01/2022 8:48:57 AM Diagnostic End Time: 28/01/2022 9:51:21 AM
Customer: UAE Service Engineer: Kanyakorn S.
Address: Contact Details: 026376363#1

Instrument Configuration

Configuration:

Serial Number: MY13160001 Turret Type: Automatic
Instrument Model: Varian AA140/240/280 Number Of Lamps: 4
Flame Instrument: True Mono Type: Automatic
Furnace Instrument: True Gasbox Type: Y Gas Box
Zeeman Present: False Auto Burner Adjuster: False
Internal Zeeman: False Mains Frequency: 50
Internal UltrAA: False Firmware Version: 2.11
Optics Type: Double Beam Photomultiplier Type: Normal(900nm)
D2 BG Correction Fitted: True PWB Version: 45
Boot Block Version: 1.09

EEPROM Data:

Instrument Run Hours: 48968.684 D2 Run Hours: 38036.500
Zero Wavelength Offset: 30.071 D2 Serial Number: not set
Mono Correction: 0.755 D2 Install Date: 01/01/1970
Flame Hours: 23619.166 D2 Original Intensity: 1.000
D2 Last Intensity: 475.000

Frequency:

Averaging Period: 30.0
Datapoint Count: 20
Upper Limit: 51.00 Highest Measured Frequency: 50.00
Average Frequency: 50.00
Lower Limit: 49.00 Lowest Measured Frequency: 50.00

Result: **Passed**

Service Completion

Service request number: 6004901032 Date service completed: 28 Jan 2022
Agilent signature: Kanyakorn S. Customer signature: Chomthanan A.
Total number of pages in this document: 13

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021



Report Generated At: 28/01/2022 9:52:26 AM

1



Power Supply:

Averaging Period: 30.0
Datapoint Count: 20



	Lower Limit (V)	Actual (V)	Upper Limit (V)	Result:
12.00 V Rail	10.80	12.20	13.20	Passed
-12.00 V Rail	-13.20	-11.90	-10.80	Passed
5.00 V Rail	4.50	5.01	5.50	Passed
310.00 V Rail	279.00	320.00	341.00	Passed

Optics

Beam Balance:

Lamp Type: Copper
Lamp Socket Used: 3

Peak Selected: 324.80
Lamp Alignment: Performed

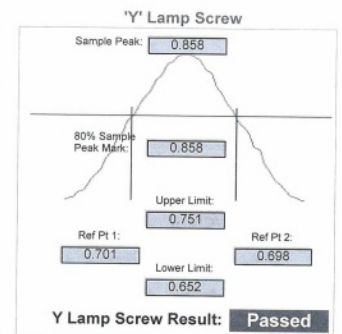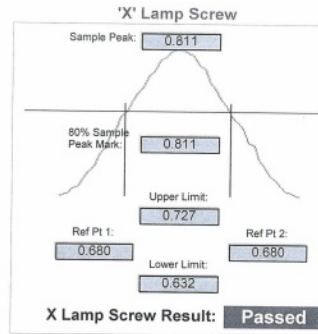

Grating Squareness:

Lamp Element(s): Copper
Lamp Turret Position: 3
Lamp Current(mA): 4.00
Slit Width(nm): 0.5
1st Order Wavelength(nm): 324.80
Lamp Alignment: Performed

	Lower Limit (nm)	Actual (nm)	Upper Limit (nm)	Result:
Zero Order	-0.10	0.00	0.10	Passed
First Order	324.45	324.75	325.15	Passed
Second Order	649.23	649.52	649.97	Passed

Wavelength Repeatability:

Lamp Used: Copper
Peak Used(nm): 324.750
Connected to Socket: 3
Lamp Current(mA): 4
Slit Width(nm): 0.2
Slit Height: Normal
Lamp Alignment: Performed

Lower Limit(nm)	Upper Limit(nm)
(Approach from Zero Order)	(Approach from end)
Sample 1: 324.823	Sample 2: 324.823
Sample 3: 324.823	Sample 4: 324.819
Sample 5: 324.819	Sample 6: 324.819
Sample 7: 324.819	Sample 8: 324.819
Sample 9: 324.819	Sample 10: 324.819

Mean: 324.821
Standard Deviation: 0.002
Result: Passed

Mechanical

Wavelength Drive:

Passed

Slit Drive:

Passed

Turret Drive:

Passed

Auto Burner Adjuster Drive:

Untested

Miscellaneous

Signal Processing Linearity:

Calculate Mode: New Calc Mode

	Lower Limit	Actual	Upper Limit	Result:
S0	114	260	297	Passed
S1	156	164	191	Passed
S2	271	295	332	Passed
S3	474	505	579	Passed
S4	825	913	1008	Passed
S5	1435	1519	1754	Passed
S6	2498	2753	3053	Passed
S7	4347	4724	5313	Passed

Interlocks:

Burner Fitted: Working	Flame Detect: Working
N2O Burner Fitted: Untested	GCU Active: Working
Flame Shield Closed: Working	Oxidant Pressure: Working
Gas Control Fitted: Untested	Oxidant Changeover: Untested
Pressure Release Bung Fitted: Working	Ignition: Working
Liquid Trap Fitted: Working	

Auto Lamp Recognition:

Lamp 1: 12 - Chromium (Cr) Lamp 5: Not Supported
Lamp 2: Uncoded Lamp/Not Connected Lamp 6: Not Supported
Lamp 3: 14 - Copper (Cu) Lamp 7: Not Supported
Lamp 4: 87 - Silver/Cadmium/Lead/Zinc(UltrAA) (Ag/C Lamp 8: Not Supported

Result: **Passed**

GTA Temperature Monitoring:

Not Performed

Notes:

Signatures:

Chonthanan A. 28 Jan 2022 Kanyakorn S. 28 Jan 2022
UAE Date Kanyakorn S. Date

Sequential by time report

28/01/2022 11:03 AM
Page 1 of 1

SpectrAA

Analyst
Date Started 28/01/2022 9:56 AM GMT: 28/01/2022 2:56 AM
Worksheet PM 28-01-2022 01-Send/ly test
Comment
Methods Cu
Computer name HEM-212
Serial Number: MY13160001

Method: Cu (Flame)

Sample ID	Exp Abs	%RSD	Mean Abs			
Cu 5 ppm	0.500	0.4	0.5004			
Readings						
	0.5042	0.5002	0.4979	0.5006	0.5021	0.4980
	0.5019	0.5004	0.4988	0.4996	28/01/2022	

Report Generated At: 28/01/2022 9:52:26 AM

6

SVD Results Report SVD
เอกสารไม่ควบคุม

Sequential by time report

28/01/2022 12:22 PM
Page 1 of 1

SpectrAA

Analyst
Date Started 28/01/2022 11:04 AM GMT: 28/01/2022 4:04 AM
Worksheet As Hydride
Comment
Methods As
Computer name HEM-212
Serial Number: MY13160001

Method: As (Vapor)

Sample ID	Exp Abs	%RSD	Mean Abs			
As 10 ppb	0.22	1.1	0.2156			
Readings						
	0.2188	0.2171	0.2150	0.2145	0.2125	28/01/2022

SPC Calibration Center

SERT
Part of DKSH Group

Certificate of Calibration

Equipment: CONDUCTIVITY METER Certificate No.: C24210091
Model: Lab955 Issued Date: 29 March 2021
Serial No. (or ID.): 16300356 Job No.: KSPR2104894
Manufacturer: SI Analytics Page: 1 of 2
Electrode Serial No. 16070067 Model: LF413T Brand: SI Analytics
Condition: In Condition

Customer: United Analyst and Engineering Consultant Company Limited
3 Soi Udomsuk 41 Sukhumvit Road,
Bangchak, Prakanong, Bangkok 10260 Thailand

Environment Condition: Temperature 23 °C ± 2 °C
Humidity 50 %RH ± 15 %RH

Calibration Place: Environment Laboratory, SPC RT Co., Ltd.
1194 Soi Wachirathamsathit 57, Sukhumvit 101/1 Rd.,
Bangchak, Prakanong, Bangkok 10260 Thailand

Calibration By: Mr.Imron Ama
Calibration Date: 29 March 2021
The Method used: In house method, SPCC-WI-49, base on ASTM D 1125-14 and D 5391-14
Traceability: This certificate is traceable to the CRM maintained by DAKKS/DKD calibration laboratory through Radiometer Analytical Co., Ltd. Certificate No. 1561, 1515, 1377

(Mr. Imron Ama)
Person in chargeบริษัท เอสพีซี อาร์ที จำกัด
SPC RT Co., Ltd.(Mr. Dumrong Boonsopon)
Authorized signatory

This certificate is issued the units of measurement according to the International System of Units (SI). It provides traceability of measurement to international or national standard or other recognized national standard laboratories.
The measurement uncertainty stated is the expanded uncertainty which is obtained from the standard uncertainty multiplied by the coverage factor ($k=2$) to provide a level of confidence of approximately 95%. It is determined in accordance with the Guide to Expression of Uncertainty in Measurement (GUM).
These results may be affected by deviations from specified conditions. The results relate only to the items tested, calibrated or sampled. The report shall not be reproduced except in full without approval of SPC RT Co., Ltd.

Calibration Results:

Before Adjustment

Standard Conductivity Solution	Unit Under Calibration Reading	Correction	Coverage Factor (k)	Uncertainty (±)
24.97 μS/cm	26.7 μS/cm	-1.73 μS/cm	2.00	0.52 μS/cm
1408.3 μS/cm	1439 μS/cm	-30.7 μS/cm	2.00	7.8 μS/cm
111.31 mS/cm	112.4 mS/cm	-1.09 mS/cm	2.00	0.58 mS/cm

After Adjustment ; at 1408.3 μS/cm

Standard Conductivity Solution	Unit Under Calibration Reading	Correction	Coverage Factor (k)	Uncertainty (±)
24.97 μS/cm	25.8 μS/cm	-0.83 μS/cm	2.00	0.52 μS/cm
1408.3 μS/cm	1410 μS/cm	-1.7 μS/cm	2.00	7.8 μS/cm
111.31 mS/cm	110.1 mS/cm	1.21 mS/cm	2.00	0.58 mS/cm

The End of Certificate

Calibration Certificate

Certificate No.: 2201793-001-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road, Bangchack, Prakhonong, Bangkok 10260

Page 1 of 5

Equipment: pH Meter
Manufacturer: METTLER TOLEDO
Model: SevenEasy pH
Serial No.: 1231155210
ID No.: UAE.WAT.010/2553
Order No.: 2201793
Operation No.: 2201793-001
Date of Receipt: 21 February 2022
Date of Calibration: 1 March 2022

Calibrated by Mr.Pheraphat Tuanjit Scientist
Approved by *P. Jenghambit* (Mr.Nuttapol Niyomchart) Specialist, Division of Calibration Laboratory
Date of Issue: 1 March 2022 Responsible for the Technical Management Team

The uncertainties are for a confidence probability of approximately 95%.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 00 Date: 14-12-61

เอกสารไม่ควบคุม

Calibration Report

Certificate No.: 2201793-001-01
Equipment: pH Meter
Resolution: 0.01 pH : 1 mV
Manufacturer: METTLER TOLEDO
Model: SevenEasy pH
Serial No.: 1231155210
Type: Bench top
ID No.: UAE.WAT.010/2553
Date of Calibration: 1 March 2022 Page 2 of 5
Location: Chemical Calibration Laboratory, NATIONAL FOOD INSTITUTE
Environment Condition: Ambient Temperature: (23.5 ± 1.5) °C Relative Humidity: (53 ± 5) %
Condition of Equipment: Good Condition
Condition of this Results of Calibration
1. Calibration Method
In house method ; W-CC-002 based on direct measurement by using standard voltage calibrator and certified reference material (CRM)
2. Reference Standards / Certified Reference Material
Instruments Serial / ID No. Manufacturer Certificate No. Due Date
2.1 DC Voltage Calibrator 2709007 Fluke SCL-21F-0687 24 June 2022
2.2 Digital Thermometer 2709007 Fluke CG-640599-01 30 October 2022
2.3 Thermo-Hygro Meter NFI.BTH004/18 PONPE QR22-0195 27 January 2023
Certified Reference Material Lot No. Manufacturer Ref N Expire Date
2.4 pH buffer 4.008 (Primary pH buffer Solution) 741339 CPAchem PH216.L5 19 April 2023
2.5 pH buffer 6.865 (Primary pH buffer Solution) 741340 CPAchem PH217.L5 19 April 2023
2.6 pH buffer 10.01 (Primary pH buffer Solution) 741342 CPAchem PH220.L5 19 April 2022
2.7 pH buffer 7.00 (Standard pH buffer Solution) 735836 CPAchem PH107.L5 16 March 2022
3. This certification is traceable to The International System of Unit (SI Unit)
3.1 Instruments No.2.1 through NSC-TISI-TIS 17025 Laboratory Accreditation of Calibration No.0075
3.2 Instruments No.2.2 through NSC-TISI-TIS 17025 Laboratory Accreditation of Calibration No.0061
3.3 Instruments No.2.3 through NSC-TISI-TIS 17025 Laboratory Accreditation of Calibration No.0292
3.4 Certified Reference Material No. 2.4 to 2.6 traceable to Primary measurement method- Hamed call using calibrated thermometer, barometer, and nanopvometer. The Standard Solution preparation and certified by CPAchem Ltd is accredited to ISO 17034 and ISO/IEC 17025
3.5 Certified Reference Material No. 2.7 traceable to BIM RefN HI-7 LotN 30.04.2020; BIM RefN HI-9 LotN 28.05.2020; BIM RefN HI-8 LotN 30.04.2020; BIM RefN HI-10 LotN 28.05.2020. The Standard Solution preparation and certified by CPAchem Ltd is accredited to ISO 17034 and ISO/IEC 17025
4. This certificate was certified only for the instrument we calibrated.
5. This result of calibration was found accurate as shown on date and place of calibration only.

P. Jenghambit
1 March 2022

เอกสารไม่ควบคุม

Calibration Report

Certificate No.: 2201793-001-01
Equipment: pH Meter
Resolution: 0.01 pH : 1 mV
Manufacturer: METTLER TOLEDO
Model: SevenEasy pH
Serial No.: 1231155210
Type: Bench top
ID No.: UAE.WAT.010/2553
Date of Calibration: 1 March 2022 Page 3 of 5
Calibration Results:
1. Calibration of pH Meter (Manual Temperature Compensation at 25 °C)
Nominal pH DC Voltage Standard (mV) Average Indicator Reading Uncertainty (±mV) Coverage Factor (k)
mV pH
0.00 414.117 414 0.00 0.58 2.00
2.00 295.811 296 2.00 0.58 2.00
4.00 177.462 178 4.00 0.58 2.00
6.00 59.159 59 6.00 0.58 2.00
7.00 -0.001 0 7.00 0.58 2.00
8.00 -59.159 -59 8.00 0.58 2.00
10.00 -177.463 -177 10.00 0.58 2.00
12.00 -295.812 -296 12.00 0.58 2.00
14.00 -414.119 -414 14.00 0.58 2.00
2. Calibration of pH Meter with Electrode (Manual Temperature Compensation at 25 °C)
Equipment: pH Electrode Type: Combined Electrode
Manufacturer: METTLER TOLEDO Model: InLabSolds
Serial No.: 1156882 ID.No. N/A
Performance of Electrode system (Three-Point Calibration at pH4, pH7 and pH10)
Certified Value @25 °C (pH) Average Indicator Reading Relative Slope (%) Uncertainty (± pH) Coverage Factor (k)
pH mV
4.008 4.00 180 96.25 0.0076 2.00
6.866 6.88 16 - 0.0079 2.00
10.012 10.01 -162 96.13 0.0094 2.00
6.865 7.00 9 - 0.0097 2.00

P. Jenghambit
1 March 2022

เอกสารไม่ควบคุม

Calibration Report

Certificate No.: 2201793-001-01
Equipment: Digital Thermometer with RTD (pH Meter)
Resolution: 0.1 °C Model: SevenEasy pH
Serial No.: 1231155210 ID No.: UAE.WAT.010/2553
Manufacturer: METTLER TOLEDO
Date of Calibration: 1 March 2022 Page 4 of 5

Location: Chemical Calibration Laboratory, NATIONAL FOOD INSTITUTE
Environment Condition: Ambient Temperature 24 °C ± 1 °C
Relative Humidity 53 % ± 2 %

Condition of this results of Calibration:

- Calibration Method :
 - In house method: W-TE-025 by comparison with standard thermometer.
 - The Calibration is determined by comparing with a known temperature from a standard resistance thermometer.
 - The temperature scale in use at this laboratory is the International Temperature scale of 1990 (ITS-90).

2. Reference Standard Instrument :

Instrument	Model	Serial No.	Certificate No.	Due Date	Through
HANDHELD THERMOMETER	1523	2118154	PSL-T 0851/64	03-Jun-22	TISTR
Platinum Resistance Thermometer (PRT)	5627A	877332			

Support Equipment : - Low Temperature Bath (ISOCAL-6), Model: Europa-6 Plus Basic, S/N: 341592/2

3. This certificate is traceable to International System of Units (SI Units).

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

6. Condition of Calibrated item :

Good

7. Result of Calibration :

☒ X

Without adjustment

☐

After adjustment

P. Pongphat
1 March 2022

F-CS-012 Revision: 00 Date: 14-12-61

เอกสารไม่ควบคุม

Calibration Report

Certificate No.: 2201793-001-01
Equipment: Digital Thermometer with RTD (pH Meter)
Resolution: 0.1 °C Model: SevenEasy pH
Serial No.: 1231155210 ID No.: UAE.WAT.010/2553
Manufacturer: METTLER TOLEDO
Date of Calibration: 1 March 2022 Page 5 of 5

Calibration point: 15.0, 25.0 and 35.0 °C
Calibration result:

- The probe was immersed in liquid bath or dry bath to a minimum depth of 100 mm.
- Description of probe, model : N/A S/N : N/A
- Dimension of probe : Diameter 4 mm., Length 100 mm.,
- Sheath material : Stainless Steel

UUC* Reading (°C)	Standard Temperature (°C)	Correction Value (°C)	Uncertainty ± (°C)
15.1	15.006	-0.1	0.099
25.1	25.004	-0.1	0.099
35.1	35.003	-0.1	0.099

P. Pongphat
1 March 2022

Note

* UUC* : Unit Under Calibration

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor $k=2$, providing a level of confidence of approximately 95 %.

----- End -----

F-CS-012 Revision: 00 Date: 14-12-61

เอกสารไม่ควบคุม

Calibration Certificate

Certificate No.: 2202097-001-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road,
Bangchack, Prakhonong, Bangkok 10260

Page 1 of 5

Equipment: pH Meter
Manufacturer: HANNA INSTRUMENTS
Model: HI 2211
Serial No.: 08165345
ID No.: UAE.WAT.004/2556
Order No.: 2202097
Operation No.: 2202097-001
Date of Receipt: 11 March 2022
Date of Calibration: 16 March 2022

Calibrated by Mr.Manas Somsak Specialist
Approved by *P. Pongphat* (Mr.Pheraphat Tuanjit)
Manager, Division of Calibration Laboratory
Date of Issue: 21 March 2022 Responsible for the Technical Management Team

The uncertainties are for a confidence probability of approximately 95%.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 00 Date: 14-12-61

เอกสารไม่ควบคุม

Calibration Report

Certificate No.: 2202097-001-01
Equipment: pH Meter
Resolution: 0.01 pH ; 0.1 mV
Manufacturer: HANNA INSTRUMENTS Model: HI 2211
Serial No.: 08165345 Type: Bench top
ID No.: UAE.WAT.004/2556

Date of Calibration: 16 March 2022 Page 2 of 5

Location: Chemical Calibration Laboratory, National Food Institute.
Environment Condition: Ambient Temperature: (23.0 ± 1.5) °C Relative Humidity: (49.5 ± 5) %
Condition of Equipment: Good Condition

Condition of this Results of Calibration

- Calibration Method : In house method : W-CC-002 based on direct measurement by using standard voltage calibrator and certified reference material (CRM)
- Reference Standards / Certified Reference Material

Instruments	Serial / ID No.	Manufacturer	Certificate No.	Due Date
2.1 DC Voltage Calibrator	2709007	Fluke	SCL-21F-0687	24 June 2022
2.2 Digital Thermometer	2709007	Fluke	CC-640599-01	30 October 2022
2.3 Thermo-Hygro Meter	ana.kh.BTH 005/58	PONPE	QR21-2787	15 November 2022

Certified Reference Material	Lot No.	Manufacturer	Ref N	Expiry Date
2.4 pH buffer 4.008 (Primary pH buffer Solution)	780012	CPAchem	PH216.L5	21 November 2023
2.5 pH buffer 6.865 (Primary pH buffer Solution)	780013	CPAchem	PH217.L5	21 November 2023
2.6 pH buffer 10.01 (Primary pH buffer Solution)	780015	CPAchem	PH220.L5	21 November 2022
2.7 pH buffer 7.00 (Standard pH buffer Solution)	776840	CPAchem	PH107.L5	8 November 2022

- This certificate is traceable to The International System of Unit (SI Unit)
 - Instruments No.2.1 through NSC-TISI-TIS 17025 Laboratory Accreditation of Calibration No.0075
 - Instruments No.2.2 through NSC-TISI-TIS 17025 Laboratory Accreditation of Calibration No.0061
 - Instruments No.2.3 through NSC-TISI-TIS 17025 Laboratory Accreditation of Calibration No.0292
 - Certified Reference Material No. 2.4 to 2.6 traceable to Primary measurement method- Harned cell using calibrated thermometer, barometer, and nanovoltmeter. The Standard Solution preparation and certified by CPAchem Ltd is accredited to ISO 17034 and ISO/IEC 17025
 - Certified Reference Material No. 2.7 traceable to BIM RefN HI-7 LotN 30.04.2020; BIM RefN HI-9 LotN 28.05.2020; BIM RefN HI-8 LotN 30.04.2020; BIM RefN HI-10 LotN 28.05.2020. The Standard Solution preparation and certified by CPAchem Ltd is accredited to ISO 17034 and ISO/IEC 17025

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

P. Pongphat

F-CS-012 Revision: 00 Date: 14-12-61

เอกสารไม่ควบคุม

Calibration Report

Certificate No.: 2202097-001-01
Equipment: pH Meter
Resolution: 0.01 pH ; 0.1 mV
Manufacturer: HANNA INSTRUMENTS Model: HI 2211
Serial No.: 08165345 Type: Bench top
ID No.: UAE.WAT.004/2556

Date of Calibration: 16 March 2022 Page 3 of 5

Calibration Results:

1. Calibration of pH Meter

(Manual Temperature Compensation at 25 °C)

Nominal pH	DC Voltage Standard (mV)	Average Indicator Reading		Uncertainty (mV)	Coverage Factor (K)
		mV	pH		
0	414.117	414	0.00	0.58	2.00
2	295.811	295.7	2.00	0.063	2.00
4	177.462	177.4	4.00	0.063	2.00
6	59.159	59.2	6.00	0.063	2.00
7	-0.001	0.1	7.00	0.063	2.00
8	-59.159	-59.1	8.00	0.063	2.00
10	-177.463	-177.3	10.00	0.063	2.00
12	-295.812	-295.6	12.00	0.063	2.00
14	-414.119	-414	14.00	0.58	2.00

2. Calibration of pH Meter with Electrode (Manual Temperature Compensation at 25 °C)

Equipment: pH Electrode Type: Combined Electrode
Manufacturer: METTLER TOLEDO Model: LE420
Serial No.: 1142802 ID No.: N/A

Performance of Electrode system

(Three-Point Calibration at pH4, pH7 and pH10)

Certified Value (25 °C (pH))	Average Indicator Reading		Relative Slope (%)	Uncertainty (± pH)	Coverage Factor (K)
	pH	mV			
4.008	4.01	180.5	99.3	0.0071	2.00
6.866	6.87	12.5	-	0.0074	2.00
10.015	10.01	-171.5	99.1	0.0090	2.00
6.983	6.98	5.2	-	0.0092	2.00

F-CS-012 Revision: 00 Date: 14-12-61

เอกสารไม่ควบคุม

Calibration Report

Certificate No.: 2202097-001-01
Equipment: Digital Thermometer with RTD (pH Meter)
Resolution: 0.1 °C Model: HI 2211
Serial No.: 08165345 ID No.: UAE.WAT.004/2556
Manufacturer: HANNA INSTRUMENTS

Date of Calibration: 16 March 2022 Page 4 of 5

Location: Chemical Calibration Laboratory, National Food Institute.

Environment Condition: Ambient Temperature (23.0 ± 1.0) °C
Relative Humidity (50 ± 4) %

Condition of this results of Calibration:

- Calibration Method : - In house method: W-TE-025 by comparison with standard thermometer.
- The Calibration is determined by comparing with a known temperature from a standard resistance thermometer.
- The temperature scale in use at this laboratory is the International Temperature scale of 1990 (ITS-90).
- Reference Standard Instrument :

Instrument	Model	Serial No.	Certificate No.	Due Date	Through
HANDHELD THERMOMETER	1523	2118154	PSL-T 8851/64	24-Jun-22	TISTR
Platinum Resistance Thermometer (PRT)	5627A	877332			

Support Equipment : - Low Temperature Bath (ISOCAL-8), Model: Europa-6 Plus Basic, S/N: 341592/2

- This certificate is traceable to International System of Units (SI Units).
- This certificate was certified only for the instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.
- Condition of Calibrated item : Good
5. Result of Calibration : ☒ Without adjustment ☐ After adjustment

F-CS-012 Revision: 00 Date: 14-12-61

เอกสารไม่ควบคุม

Calibration Report

Certificate No.: 2202097-001-01
Equipment: Digital Thermometer with RTD (pH Meter)
Resolution: 0.1 °C Model: HI 2211
Serial No.: 08165345 ID No.: UAE.WAT.004/2556
Manufacturer: HANNA INSTRUMENTS

Date of Calibration: 16 March 2022 Page 5 of 5

Calibration point: 15.0, 25.0 and 35.0 °C

Calibration result:

- The probe was immersed in liquid bath or dry bath to a minimum depth of 100 mm.

- Description of probe, model : N/A SIN : N/A

Dimension of probe : Diameter 3.5 mm., Length 100 mm.,

Sheath material : Stainless Steel

UUC* Reading (°C)	Standard Temperature (°C)	Correction Value (°C)	Uncertainty ± (°C)
15.0	15.001	0.0	0.099
25.0	25.002	0.0	0.099
35.0	35.002	0.0	0.099

Note : - UUC* : Unit Under Calibration

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k = 2, providing a level of confidence of approximately 95 %.

----- End -----

F-CS-012 Revision: 00 Date: 14-12-61

เอกสารไม่ควบคุม

Calibration Certificate

Certificate No.: 2203120-001-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO., LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road, Bangchack, Prakhonong, Bangkok 10260

Page 1 of 3

Equipment: Electronic Balance

Manufacturer: METTLER TOLEDO

Model: AB204-S/FACT

Serial No.: 1129361010

ID No.: UAE.WAS.002/2552

Order No.: 2203120

Operation No.: 2203120-001

Date of Receipt: 1 June 2022

Date of Calibration: 1 June 2022

Calibrated by Mr.Taveesak Sellee
Scientist

Approved by (Mr.Pheraphat Tuanjit)
Manager, Division of Calibration Laboratory
Responsible for the Technical Management Team

The uncertainties are for a confidence probability of approximately 95%

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65

เอกสารไม่ควบคุม



REPORT OF CALIBRATION

Certificate No. : SP22-016

Page 3 of 5

Calibration Results : Without adjustment

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor k
420	0.0000	0.0000	0.0000	0.0028	2.00
	0.5787	0.5755	0.0032	0.0031	2.00
	1.0490	1.0436	0.0054	0.0029	2.00
	2.1900	2.1847	0.0053	0.0075	2.00
440	0.0000	0.0000	0.0000	0.0028	2.00
	0.5607	0.5588	0.0019	0.0034	2.00
	1.0247	1.0232	0.0015	0.0035	2.00
	2.1229	2.1211	0.0018	0.0082	2.00
465	0.0000	0.0000	0.0000	0.0028	2.00
	0.5236	0.5197	0.0039	0.0029	2.00
	0.9634	0.9625	0.0009	0.0028	2.00
	1.9763	1.9752	0.0011	0.0070	2.00
546.1	0.0000	-0.0001	0.0001	0.0028	2.00
	0.5191	0.5171	0.0020	0.0031	2.00
	1.0003	0.9984	0.0019	0.0033	2.00
	1.9987	1.9946	0.0041	0.0084	2.00
590	0.0000	0.0000	0.0000	0.0028	2.00
	0.5523	0.5509	0.0014	0.0030	2.00
	1.0809	1.0799	0.0010	0.0029	2.00
	2.0391	2.0329	0.0062	0.0080	2.00
635	0.0000	0.0000	0.0000	0.0028	2.00
	0.5601	0.5584	0.0017	0.0031	2.00
	1.0512	1.0498	0.0014	0.0029	2.00
	1.9294	1.9265	0.0029	0.0082	2.00

เอกสารไม่ควบคุม

FM-708-02 R01 1/11/2021



REPORT OF CALIBRATION

Certificate No. : SP22-016

Page 4 of 5

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor k
235	0.0000	0.0001	-0.0001	0.0050	2.00
	0.7478	0.7421	0.0057	0.0056	2.00
257	0.0000	0.0000	0.0000	0.0050	2.00
	0.8686	0.8619	0.0067	0.0059	2.00
313	0.0000	0.0000	0.0000	0.0050	2.00
	0.2912	0.2896	0.0016	0.0051	2.00
350	0.0000	0.0000	0.0000	0.0050	2.00
	0.6448	0.6403	0.0045	0.0055	2.00

เอกสารไม่ควบคุม

FM-708-02 R01 1/11/2021



REPORT OF CALIBRATION

Certificate No. : SP22-016

Page 5 of 5

Wavelength Accuracy :

CRMs Values (nm.)	UUC Reading (nm.)	Correction (nm.)	Uncertainty (nm.)	Coverage factor k
241.72	242.0	-0.28	0.18	2.00
279.45	279.5	-0.05	0.18	2.00
287.81	287.5	0.31	0.18	2.00
334.06	333.5	0.56	0.18	2.00
360.93	360.5	0.43	0.18	2.00
418.59	418.0	0.59	0.18	2.00
445.94	445.4	0.54	0.18	2.00
453.66	453.2	0.46	0.18	2.00
460.02	459.7	0.32	0.18	2.00
536.59	536.2	0.39	0.18	2.00
637.98	638.3	-0.32	0.18	2.00
431.38	431.0	0.38	0.18	2.00
472.50	472.5	0.00	0.18	2.00
513.47	513.5	-0.03	0.18	2.00
528.88	528.5	0.38	0.18	2.00
573.17	573.0	0.17	0.18	2.00
585.35	585.0	0.35	0.20	2.00
684.40	684.7	-0.30	0.18	2.00
740.72	740.8	-0.08	0.20	2.00
748.55	748.5	0.05	0.18	2.00
807.03	807.3	-0.27	0.18	2.00
879.28	879.0	0.28	0.18	2.00

Remark : - UUC = Unit Under Calibration

- N/A = Not Available

- The result expanded uncertainty of measurement U is stated as the standard uncertainty of measurement multiplied by the coverage factor k,

which for a normal distribution corresponds to a coverage probability of approximately 95%

- * Indicates non TISI accredited

- End of Certificate -

เอกสารไม่ควบคุม

FM-708-02 R01 1/11/2021



CERTIFICATE OF CALIBRATION

Certificate No. : SP22-007

Page 1 of 5

Customer : United Analyst and Engineering Consultant Co.,Ltd. (Head Office)

Address : 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong,

Bangkok 10260

Location of calibration : Laboratory 315

Equipment : UV-Vis Spectrophotometer

Manufacturer : Hitachi

Model : U-1900

Serial No. : 2021-064

ID No. : UAE.WAS.006/2552

Received Date : 20 January 2022

Calibration Date : 20 January 2022

Issue Date : 24 January 2022

Condition Instrument : Good

Calibrated by :

(Mr. Tanawut Rittidach)

Technical Manager

Approved by :

(Ms. Chonticha Sangngern)

Quality Manager


The calibration result is applied only to the above calibrated item and was found accurate as shown on date and place of calibration only.

The measurement capability of the laboratory and its traceability to recognized national standards and to the unit of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the DQE Services Co., Ltd.

FM-708-02 R01 1/11/2021

เอกสารไม่ควบคุม

DQE Services Co.,Ltd.
32 Soi Ladprao-Wanghin 55, Ladprao-Wanghin Rd., Ladprao, Bangkok 10230
Phone : +66 (0)2 538 2054, Email : dqeservicesinfo@gmail.com



REPORT OF CALIBRATION

Certificate No. : SP22-007 Page 2 of 5

Environment Condition : Ambient Temperature 25 ± 5 °C

Relative humidity 55 ± 20 %RH

Calibration method : In-house method CP-01 Based on ASTM E275-08

Certified Reference Materials :

Material	Serial No.	Certificate No.	Due date
Absorbance Standard set	25760	95935	22 October 2023
Absorbance Standard set	25757	95929	22 October 2023
Wavelength Standard set	25806	95916	22 October 2023
Wavelength Standard set	25758	95915	22 October 2023

Traceability This certification is traceable to the International System of Unit maintained at National -
Institute of Standards and Technology (NIST) through Starna Scientific Limited

Spectral Band Width of UUC : 4.0 nm.

Scan Speed of UUC : 200 nm/min

Scan Interval of UUC : 0.1 nm.


Resolution of UUC : Photometric 0.001 Abs.

Wavelength 0.1 nm.

FM-708-02 R01 1/11/2021

เอกสารไม่ควบคุม

DQE Services Co.,Ltd.
32 Soi Ladprao-Wanghin 55, Ladprao-Wanghin Rd., Ladprao, Bangkok 10230
Phone : +66 (0)2 538 2054, Email : dqeservicesinfo@gmail.com



REPORT OF CALIBRATION

Certificate No. : SP22-007 Page 3 of 5

Calibration Results : Without adjustment


Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor k
420	0.0000	0.000	0.0000	0.0028	2.00
	0.5787	0.577	0.0017	0.0031	2.00
	1.0490	1.050	-0.0010	0.0029	2.00
	2.1900	2.183	0.0070	0.0080	2.00
440	0.0000	0.000	0.0000	0.0028	2.00
	0.5607	0.560	0.0007	0.0034	2.00
	1.0247	1.023	0.0017	0.0035	2.00
	2.1229	2.118	0.0049	0.0079	2.00
465	0.0000	0.000	0.0000	0.0028	2.00
	0.5236	0.521	0.0026	0.0030	2.00
	0.9634	0.963	0.0004	0.0029	2.00
	1.9763	1.974	0.0023	0.0070	2.00
546.1	0.0000	0.000	0.0000	0.0028	2.00
	0.5191	0.518	0.0011	0.0031	2.00
	1.0003	1.000	0.0003	0.0033	2.00
	1.9987	1.996	0.0027	0.0084	2.00
590	0.0000	0.000	0.0000	0.0028	2.00
	0.5523	0.552	0.0003	0.0030	2.00
	1.0809	1.082	-0.0011	0.0030	2.00
	2.0391	2.033	0.0061	0.0079	2.00
635	0.0000	0.000	0.0000	0.0028	2.00
	0.5601	0.562	-0.0019	0.0031	2.00
	1.0512	1.052	-0.0008	0.0030	2.00
	1.9294	1.925	0.0044	0.0079	2.00

FM-708-02 R01 1/11/2021

เอกสารไม่ควบคุม

DQE Services Co.,Ltd.
32 Soi Ladprao-Wanghin 55, Ladprao-Wanghin Rd., Ladprao, Bangkok 10230
Phone : +66 (0)2 538 2054, Email : dqeservicesinfo@gmail.com



REPORT OF CALIBRATION

Certificate No. : SP22-007 Page 4 of 5


Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor k
235	0.0000	0.000	0.0000	0.0050	2.00
	0.7478	0.746	0.0018	0.0057	2.00
257	0.0000	0.000	0.0000	0.0050	2.00
	0.8686	0.861	0.0076	0.0059	2.00
313	0.0000	0.000	0.0000	0.0050	2.00
	0.2912	0.291	0.0002	0.0051	2.00
350	0.0000	0.000	0.0000	0.0050	2.00
	0.6448	0.638	0.0068	0.0055	2.00

FM-708-02 R01 1/11/2021

เอกสารไม่ควบคุม

DQE Services Co.,Ltd.
32 Soi Ladprao-Wanghin 55, Ladprao-Wanghin Rd., Ladprao, Bangkok 10230
Phone : +66 (0)2 538 2054, Email : dqeservicesinfo@gmail.com



REPORT OF CALIBRATION

Certificate No. : SP22-007 Page 5 of 5

Wavelength Accuracy :

CRMs Values (nm.)	UUC Reading (nm.)	Correction (nm.)	Uncertainty (nm.)	Coverage factor k
241.54	240.8	0.74	0.18	2.00
279.40	278.5	0.90	0.18	2.00
288.70	288.0	0.70	0.18	2.00
334.22	333.5	0.72	0.18	2.00
361.26	360.5	0.76	0.18	2.00
418.48	418.0	0.48	0.18	2.00
446.70	446.0	0.70	0.18	2.00
453.20	453.0	0.20	0.18	2.00
460.06	459.5	0.56	0.18	2.00
536.90	536.0	0.90	0.18	2.00
637.94	637.2	0.74	0.18	2.00
440.74	440.0	0.74	0.18	2.00
472.22	471.6	0.62	0.18	2.00
513.70	513.0	0.70	0.18	2.00
528.72	528.0	0.72	0.18	2.00
574.60	573.8	0.80	0.18	2.00
585.48	584.6	0.88	0.20	2.00
684.63	684.0	0.63	0.18	2.00
740.27	739.8	0.47	0.20	2.00
748.28	747.8	0.48	0.18	2.00
807.16	806.4	0.76	0.18	2.00
879.70	878.8	0.90	0.18	2.00

Remark : - UUC = Unit Under Calibration

- N/A = Not Available

- The result expanded uncertainty of measurement U is stated as the standard uncertainty of measurement multiplied by the coverage factor k,
which for a normal distribution corresponds to a coverage probability of approximately 95%

- * Indicates non TISI accredited

- End of Certificate -

FM-708-02 R01 1/11/2021

เอกสารไม่ควบคุม

DQE Services Co.,Ltd.
32 Soi Ladprao-Wanghin 55, Ladprao-Wanghin Rd., Ladprao, Bangkok 10230
Phone : +66 (0)2 538 2054, Email : dqeservicesinfo@gmail.com

CERTIFICATE OF CALIBRATION

Certificate No. : SP22-008 Page 1 of 5

Customer : United Analyst and Engineering Consultant Co.,Ltd. (Head Office)

Address : 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong,
Bangkok 10260

Location of calibration : Laboratory 213

Equipment : UV-Vis Spectrophotometer

Manufacturer : Hitachi

Model : U-2900

Serial No. : 21E22-009



ID No. : UAE.WAT.051/2564

Received Date : 20 January 2022

Calibration Date : 20 January 2022

Issue Date : 24 January 2022

Condition Instrument : Good

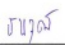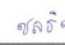Calibrated by :  Approved by : 
(Mr. Tanawut Rittidach) (Ms. Chonthicha Sangngern)
Technical Manager Quality Manager

The calibration result is applied only to the above calibrated item and was found accurate as shown on date and place of calibration only.

The measurement capability of the laboratory and its traceability to recognized national standards and to the unit of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the DQE Services Co., Ltd.

FM-708-02 R01 1/11/2021

เอกสารไม่ควบคุม

DQE Services Co.,Ltd.
32 Soi Ladprao-Wanghin 55, Ladprao-Wanghin Rd., Ladprao, Bangkok 10230
Phone : +66 (0)2 538 2054, Email : dqeservicesinfo@gmail.com

REPORT OF CALIBRATION

Certificate No. : SP22-008 Page 2 of 5

Environment Condition : Ambient Temperature 25 ± 5 °C
Relative humidity 55 ± 20 %RH

Calibration method : In-house method CP-01 Based on ASTM E275-08

Certified Reference Materials :

Material	Serial No.	Certificate No.	Due date
Absorbance Standard set	25760	95935	22 October 2023
Absorbance Standard set	25757	95929	22 October 2023
Wavelength Standard set	25806	95916	22 October 2023
Wavelength Standard set	25758	95915	22 October 2023

Traceability This certification is traceable to the International System of Unit maintained at National -
Institute of Standards and Technology (NIST) through Starna Scientific Limited

Spectral Band Width of UUC : 1.5 nm.

Scan Speed of UUC : 200 nm/min

Scan Interval of UUC : 0.1 nm.

Resolution of UUC : Photometric 0.001 Abs.
Wavelength 0.1 nm.

FM-708-02 R01 1/11/2021

เอกสารไม่ควบคุม

DQE Services Co.,Ltd.
32 Soi Ladprao-Wanghin 55, Ladprao-Wanghin Rd., Ladprao, Bangkok 10230
Phone : +66 (0)2 538 2054, Email : dqeservicesinfo@gmail.com

REPORT OF CALIBRATION

Certificate No. : SP22-008 Page 3 of 5

Calibration Results : Without adjustment

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor k
420	0.0000	0.000	0.0000	0.0028	2.00
	0.5787	0.576	0.0027	0.0031	2.00
	1.0490	1.046	0.0030	0.0029	2.00
	2.1900	2.182	0.0080	0.0075	2.00
440	0.0000	0.000	0.0000	0.0028	2.00
	0.5607	0.559	0.0017	0.0034	2.00
	1.0247	1.023	0.0017	0.0035	2.00
	2.1229	2.116	0.0069	0.0079	2.00
465	0.0000	0.000	0.0000	0.0028	2.00
	0.5236	0.521	0.0026	0.0030	2.00
	0.9634	0.962	0.0014	0.0029	2.00
	1.9763	1.970	0.0063	0.0070	2.00
546.1	0.0000	0.000	0.0000	0.0028	2.00
	0.5191	0.519	0.0001	0.0031	2.00
	1.0003	0.999	0.0013	0.0033	2.00
	1.9987	1.992	0.0067	0.0084	2.00
590	0.0000	0.000	0.0000	0.0028	2.00
	0.5523	0.552	0.0003	0.0030	2.00
	1.0809	1.080	0.0009	0.0030	2.00
	2.0391	2.031	0.0081	0.0079	2.00
635	0.0000	0.000	0.0000	0.0028	2.00
	0.5601	0.560	0.0001	0.0031	2.00
	1.0512	1.052	-0.0008	0.0030	2.00
	1.9294	1.922	0.0074	0.0079	2.00

FM-708-02 R01 1/11/2021

เอกสารไม่ควบคุม

DQE Services Co.,Ltd.
32 Soi Ladprao-Wanghin 55, Ladprao-Wanghin Rd., Ladprao, Bangkok 10230
Phone : +66 (0)2 538 2054, Email : dqeservicesinfo@gmail.com

REPORT OF CALIBRATION

Certificate No. : SP22-008 Page 4 of 5

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor k
235	0.0000	0.000	0.0000	0.0050	2.00
	0.7478	0.747	0.0008	0.0057	2.00
257	0.0000	0.000	0.0000	0.0050	2.00
	0.8686	0.865	0.0036	0.0059	2.00
313	0.0000	0.000	0.0000	0.0050	2.00
	0.2912	0.290	0.0012	0.0051	2.00
350	0.0000	0.000	0.0000	0.0050	2.00
	0.6448	0.640	0.0048	0.0055	2.00

FM-708-02 R01 1/11/2021

เอกสารไม่ควบคุม

REPORT OF CALIBRATION

Certificate No. : SP22-008

Page 5 of 5

Wavelength Accuracy :

CRMs Values (nm.)	UUC Reading (nm.)	Correction (nm.)	Uncertainty (nm.)	Coverage factor k
241.72	241.0	0.72	0.18	2.00
279.45	279.0	0.45	0.18	2.00
287.81	287.0	0.81	0.18	2.00
334.06	333.5	0.56	0.18	2.00
360.93	360.0	0.93	0.18	2.00
418.59	418.0	0.59	0.18	2.00
445.94	445.5	0.44	0.18	2.00
453.66	453.0	0.66	0.18	2.00
460.02	459.5	0.52	0.18	2.00
536.59	536.0	0.59	0.18	2.00
637.98	637.5	0.48	0.18	2.00
431.38	431.0	0.38	0.18	2.00
472.50	472.0	0.50	0.18	2.00
513.47	513.0	0.47	0.18	2.00
528.88	528.5	0.38	0.18	2.00
573.17	573.0	0.17	0.18	2.00
585.35	585.0	0.35	0.20	2.00
684.40	684.0	0.40	0.18	2.00
740.72	740.5	0.22	0.20	2.00
748.55	748.5	0.05	0.18	2.00
807.03	807.0	0.03	0.18	2.00
879.28	879.5	-0.22	0.18	2.00

Remark : - UUC = Unit Under Calibration

- N/A = Not Available

- The result expanded uncertainty of measurement U is stated as the standard uncertainty of measurement multiplied by the coverage factor k ,

which for a normal distribution corresponds to a coverage probability of approximately 95%

- * Indicates non TISI accredited

- End of Certificate -

FM-708-02 R01 1/11/2021

เอกสารไม่ควบคุม