

ภาคผนวก 12

เอกสารการสอบเทียบเครื่องมือและ  
อุปกรณ์ตรวจวัดคุณภาพสิ่งแวดล้อม

---



## List of Instruments Certification for Water Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration	Remark
Water									
1	pH Meter	pH	Horiba	LAQUA-PH210 HA9M0047	Technology Promotion Association (Thailand-Japan)	22CH650	13 May 22	12 May 23	-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES  
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250  
TEL. 0-2717-3000-27 FAX. 0-2719-9484



Cert.No.: 22CH650

Page.: 1 of 3

## Certificate of Calibration

**Equipment :** pH Meter  
**Manufacturer :** Horiba  
**Model :** LAQUA-PH210  
**Serial No. :** HA9M0047  
**ID No. :** UAE.EFM.005/2563(EFM.pH.05/63)  
**Condition As-Received:** Used Item  
**Received Date :** 29 April 2022  
**Calibration Date :** 13 May 2022  
**Reference :** 2204-0781WSC-1  
**Submitted by :** United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road, Bangchak,  
Phrakhanong, Bangkok 10260  
**Ambient Temperature :** (25  $\pm$  2.5) °C  
**Relative Humidity :** (50  $\pm$  15) %  
**Calibration Procedure :** In - house method :  
- CP-CH5 by direct measurement with standard  
voltage calibrator and direct measurement with  
certified reference material (CRM)  
- CP-CH8 by comparison with standard thermometer

**Calibrated by :** Uthen Kankawi

**Approved by :**

(☒) Malee Butkruea  
( ) Saithip Meangmai  
( ) Warakorn Lernagatrakul

**Issue Date :** 19 May 2022

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม

A 0041259



Cert.No.: 22CH650

Page.: 2 of 3

**Condition of this calibration result**

1. Reference Standard Instrument : -

<u>Instrument</u>	<u>Serial No.</u>	<u>ID No.</u>	<u>Cert. No.</u>	<u>Due Date</u>
1) Document Process Calibrator	46530031	130RC098	21E3245	07 Oct 2022
2) Ref. Standard Thermometer	2188080	130RC044	21I1273	21 Nov 2022

This certification is traceable to the International System of Unit maintained at:-

- Traceable to National Institute of Metrology (Thailand), NIMT

2. Certified Reference Materials : The measurement results are traceable to SI through CPA chem Ltd.,  
ANSI-ASQ National Accreditation Board, Accredited No. AR-1835

<u>Buffer Solution</u>	<u>Manufacturer</u>	<u>Lot No.</u>	<u>Exp. date</u>
pH 4.008	CPA chem	788995	01 Jan 2024
pH 6.983	CPA chem	766822	04 Sep 2022
pH 10.015	CPA chem	794124	14 Feb 2023

3. This certificate is valid only to the item calibrated on date and place of calibration.

**Calibration Results****Function : mV Measurement**

Performing standard curve by Fluke at pH (4,7)(7,10)

Unit Under Calibration	Nominal Value	Standard Voltage Input	Actual Reading		Uncertainty of Measurement ( ±mV )	Coverage factor <i>k</i>
	pH	mV	mV	pH		
pH Meter S/N.: HA9M0047	4.00	177.48	177.7	4.00	0.058	2.00
	7.00	0.00	0.3	7.01	0.058	2.00
	7.00	0.00	0.3	7.01	0.058	2.00
	10.00	-177.48	-177.0	10.00	0.058	2.00

Malu

เอกสารไม่ควบคุม

a 1108698



Cert.No.: 22CH650

Page.: 3 of 3

### Calibration Results

Function : pH Measurement

Performing three buffers standard curve by using buffer nominal pH (4,7)(7,10)

Unit Under Calibration	Standard pH Buffer Solution	Actual pH Reading	Actual mV Reading ( mV )	Uncertainty of pH measurement ( $\pm$ )	Coverage factor $k$
pH Electrode S/N.: 990C0199	4.008	4.01	131.3	0.011	2.13
	6.983	6.98	-43.3	0.012	2.07
	6.983	6.99	-43.7	0.011	2.00
	10.015	10.02	-215.3	0.013	2.13

Function : Temperature Measurement

( \* ) Without adjustment

This equipment was connected with Temperature Probe;

- Model : 9652  
- Serial No. : 990C0199

Dimension of probe;

- Length : 95 mm.  
- Diameter : 12 mm.  
- Immersion Depth : 100 mm.

Calibration Point ( $^{\circ}\text{C}$ )	Standard Temperature ( $^{\circ}\text{C}$ )	UUC* Reading ( $^{\circ}\text{C}$ )	Error ( $^{\circ}\text{C}$ )	Uncertainty of measurement ( $\pm$ $^{\circ}\text{C}$ )	Coverage factor $k$
25.0	25.004	25.0	-0.004	0.13	2.00
30.0	30.005	30.0	-0.005	0.13	2.00
35.0	35.012	35.0	-0.012	0.13	2.00

Remark : - UUC\* = Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor  $k$ , providing a level of confidence of approximately 95 %.

-o0o-

Maku.

เอกสารไม่ควบคุม

a 1108697

List Certificate Instruments of Laboratory for Environmental Quality Analysis.

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*	Remark
Laboratory Instrument/Equipments (คุณภาพน้ำ)									
1	pH Meter	pH	Mettler-Toledo	Seven Easy S20 / 1230525212	National Food Institute, Ministry of Industry, Thailand	2202093-001-01	16 Mar 22	15 Mar 23	-
2	pH Meter		Mettler-Toledo	SevenCompact S220/ C113432421	National Food Institute, Ministry of Industry, Thailand	2203527-001-01	5 Jul 22	4 Jul 23	-
3	Analytical Balance (Readability 0.01 mg)	Total Dissolved Solids (TDS) Suspended Solids (SS)	Mettler-Toledo	XSR205DU / C009071872	Technology Promotion Association (Thailand-Japan)	22MM210	26 Apr 22	25 Apr 23	-
4	Hot Air Oven		Memmert	UF55 / B216.1666	Technology Promotion Association (Thailand-Japan)	22TM1490	19 Oct 22	18 Oct 23	-
5	Analytical Balance (Readability 0.1 mg)	Fat Oil and Grease (FOG)	Mettler-Toledo	XSR204 / C117635043	National Food Institute, Ministry of Industry, Thailand	2202934-001-01	13 May 22	12 May 23	-
6	BOD Incubator	Biochemical Oxygen Demand (BOD)	Arco	UC4-1320 / (UAE.WAO.015/2561)	Technology Promotion Association (Thailand-Japan)	22TM90	17 Feb 22	16 Feb 23	-
7	BOD Incubator		Arco	UR-1320 / (UAE.WAO.018/2551)	Technology Promotion Association (Thailand-Japan)	22TM305	7 Apr 22	6 Apr 23	-
8	COD Reactor (Heating Block)	Chemical Oxygen Demand (COD)	Hanna	HI839800-02 / H018500I	Hanna Instruments (Thailand) Ltd.	HIT-2209-0184	1 Mar 22	1 Mar 23	-
9	Digester Unit	Total Kjeldahl Nitrogen (TKN)	FOSS TECATOR	2520auto / 91794469	National Food Institute, Ministry of Industry, Thailand	2202361-001-01	4 Apr 22	3 Apr 23	-
10	Distillation Unit (Kjeldahl Method)		FOSS TECATOR	KT8100/ 91889052	FOSS South East Asia	6623	25 Jul 22	24 Jul 23	-

### List Certificate Instruments of Laboratory for Environmental Quality Analysis.

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*	Remark
Laboratory Instrument/Equipments (คุณภาพน้ำ)									
11	UV-VIS Spectrophotometer	Nitrate-Nitrogen Nitrite-Nitrogen	Agilent Technologies	Cary60 G6860A / MY15410009	DQE Services Co.,Ltd.	SP22-016	31 May 22	30 May 23	-
12	UV-VIS Spectrophotometer	Ammonia-Nitrogen Total Phosphorus	Hitachi	U-1900 / 2021-064	DQE Services Co.,Ltd.	SP22-007	20 Jan 22	19 Jan 23	-
13	Incubator	Total Coliform Bacteria Fecal Coliform Bacteria	Memmert	IPP 260 / V615.0187	Technology Promotion Association (Thailand-Japan)	22TM563	7 Apr 22	6 Apr 23	-
14	Incubator	E. coliform Bacteria	Memmert	IPP 260 / V618.0033	Technology Promotion Association (Thailand-Japan)	22TM503	3 May 22	2 May 23	-
15	Water Bath		Memmert	WNE 14 / L414.1407	Technology Promotion Association (Thailand-Japan)	22TM565	7 Apr 22	6 Apr 23	-
16	Water Bath		Memmert	WNE 14 / L414.1410	Technology Promotion Association (Thailand-Japan)	22TM564	7 Apr 22	6 Apr 23	-
17	Analytical Balance		Mettler-Toledo	MS603S / B007010311	Mettler-Toledo (Thailand) Ltd.	2058-096-040722-ACC-	7 Apr 22	6 Apr 23	-
18	Auto Clave		ALP	CL-40L / 802664	Technology Promotion Association (Thailand-Japan)	22TM89	17 Feb 22	16 Feb 23	-

Due Date of Calibration\* : Based on the annual calibration plan. At least 1 time per year.

## Calibration Certificate

**Certificate No.:** 2202093-001-01  
**Client name:** UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.  
**Address:** 3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchack, Prakhonong, Bangkok 10260

Page 1 of 5

**Equipment:** pH Meter  
**Manufacturer:** METTLER TOLEDO  
**Model:** SevenEasy pH  
**Serial No.:** 1230525212  
**ID No.:** UAE.WAS.003/2553  
**Order No.:** 2202093  
**Operation No.:** 2202093-001  
**Date of Receipt:** 11 March 2022  
**Date of Calibration:** 16 March 2022

**Calibrated by** Mr.Manas Somsak  
Specialist

**Approved by** (Mr. [Redacted])  
Manager, Division of Calibration Laboratory  
Responsible for the Technical Management Team

**Date of Issue:** 21 March 2022

The uncertainties are for a confidence probability of approximately 95%.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

## Calibration Report

**Certificate No.:** 2202093-001-01  
**Equipment:** pH Meter  
**Manufacturer:** METTLER TOLEDO  
**Serial No.:** 1230525212  
**ID No.:** UAE.WAS.003/2553  
**Resolution:** 0.01 pH ; 1 mV  
**Model:** SevenEasy pH  
**Type:** Bench top

**Date of Calibration:** 16 March 2022 Page 2 of 5

**Location:** Chemical Calibration Laboratory, National Food Institute.  
**Environment Condition:** Ambient Temperature: { 23.0 ± 1.5 } °C Relative Humidity: { 49.5 ± 5 } %  
**Condition of Equipment:** Good Condition  
**Condition of this Results of Calibration**

1. Calibration Method In house method : W-CC-002 based on direct measurement by using standard voltage calibrator and certified reference material (CRM)

2. Reference Standards / Certified Reference Material

Instruments	Serial / ID No.	Manufacturer	Certificate No.	Due Date
2.1 DC Voltage Calibrator	2709007	Fluke	SCL-21F-0687	24 June 2022
2.2 Digital Thermometer	2709007	Fluke	CC-640599-01	30 October 2022
2.3 Thermo-Hygro Meter	ana.ksl.BTH 005/58	PONPE	QR21-2787	15 November 2022

Certified Reference Material	Lot No.	Manufacturer	Ref N	Expire Date
2.4 pH buffer 4.008 (Primary pH buffer Solution)	780012	CPAchem	PH216.L5	21 November 2023
2.5 pH buffer 6.865 (Primary pH buffer Solution)	780013	CPAchem	PH217.L5	21 November 2023
2.6 pH buffer 10.01 (Primary pH buffer Solution)	780015	CPAchem	PH220.L5	21 November 2022
2.7 pH buffer 7.00 (Standard pH buffer Solution)	776840	CPAchem	PH107.L5	8 November 2022

3. This certification is traceable to The International System of Unit (SI Unit)

3.1 Instruments No.2.1	through	NSC-TIS-115 17025 Laboratory Accreditation of Calibration No.0075
3.2 Instruments No.2.2	through	NSC-TIS-115 17025 Laboratory Accreditation of Calibration No.0061
3.3 Instruments No.2.3	through	NSC-TIS-115 17025 Laboratory Accreditation of Calibration No.0292
3.4 Certified Reference Material No. 2.4 to 2.6	traceable to	Primary measurement method- Harned cell using calibrated thermometer, barometer, and nanovoltmeter. The Standard Solution preparation and certified by CPAchem Ltd is accredited to ISO 17034 and ISO/IEC 17025
3.5 Certified Reference Material No. 2.7	traceable to	BIM RefN HI-7 LotN 30.04.2020; BIM RefN HI-9 LotN 28.05.2020; BIM RefN HI-8 LotN 30.04.2020; BIM RefN HI-10 LotN 28.05.2020. The Standard Solution preparation and certified by CPAchem Ltd is accredited to ISO 17034 and ISO/IEC 17025

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

## Calibration Report

Certificate No.: 2202093-001-01

Equipment:

pH Meter

Resolution: 0.01 pH ; 1 mV

Manufacturer: METTLER TOLEDO

Model: SevenEasy pH

Serial No.: 1230525212

Type: Bench top

ID No.: UAE.WAS.003/2553

Date of Calibration: 16 March 2022

Page 3 of 5

Calibration Results:

1. Calibration of pH Meter ( Manual Temperature Compensation at 25 °C )

Nominal pH	DC Voltage Standard ( mV )	Average Indicator Reading		Uncertainty ( ±mV )	Coverage Factor ( k )
		mV	pH		
0	414.117	414	0.00	0.58	2.00
2	295.811	296	2.00	0.58	2.00
4	177.482	178	4.00	0.58	2.00
6	59.159	59	6.00	0.58	2.00
7	-0.001	0	7.00	0.58	2.00
8	-59.159	-59	8.00	0.58	2.00
10	-177.463	-177	10.00	0.58	2.00
12	-295.812	-296	12.00	0.58	2.00
14	-414.119	-414	14.00	0.58	2.00

2. Calibration of pH Meter with Electrode ( Manual Temperature Compensation at 25 °C )

Equipment: pH Electrode

Type: Combined Electrode

Manufacturer: METTLER TOLEDO

Model: InLab Solids

Serial No.: 9453943

ID.No. N/A

Performance of Electrode system (Three-Point Calibration at pH4, pH7 and pH10)

Certified Value @25 °C (pH)	Average Indicator Reading		Relative Slope (%)	Uncertainty ( ± pH )	Coverage Factor ( k )
	pH	mV			
4.006	4.01	172	98.1	0.0071	2.00
6.866	6.87	6	-	0.0074	2.00
10.015	10.01	-175	97.4	0.0090	2.00
6.983	6.98	-3	-	0.0092	2.00

## Calibration Report

Certificate No.: 2202093-001-01

Equipment:

Digital Thermometer with RTD (pH Meter)

Resolution: 0.1 °C

Model: SevenEasy pH

Serial No.: 1230525212

ID No.: UAE.WAS.003/2553

Manufacturer: METTLER TOLEDO

Date of Calibration: 16 March 2022

Page 4 of 5

Location:

Chemical Calibration Laboratory, National Food Institute.

Environment Condition:

Ambient Temperature ( 23.0 ± 1.0 ) °C

Relative Humidity ( 50 ± 4 ) %

Condition of this results of Calibration:

- Calibration Method :
  - In house method: W-TE-025 by comparison with standard thermometer.
  - The Calibration is determined by comparing with a known temperature from a standard resistance thermometer.
  - The temperature scale in use at this laboratory is the International Temperature scale of 1990 ( ITS-90 ).

2. Reference Standard Instrument :

Instrument	Model	Serial No.	Certificate No.	Due Date	Through
HANDHELD THERMOMETER	1523	2118154	PSL-T 0851/64	24-Jun-22	TISTR
Platinum Resistance Thermometer (PRT)	5627A	877332			

Support Equipment : - Low Temperature Bath (ISOCAL-6), Model: Europa-6 Plus Basic, S/N: 341592/2

- This certificate is traceable to International System of Units (SI Units).
- This certificate was certified only for the instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.

6. Condition of Calibrated item : Good

7. Result of Calibration :

☒ X

Without adjustment

☐

After adjustment



National Food Institute, Ministry of Industry, Thailand

2008 Soi 36, Arun Amarin Road, Bang Yi Khan Subdistrict, Bang Phlat District, Bangkok 10700, Thailand  
Tel : +66 (0) 2422 8588 Fax : +66 (0) 2422 8558 Website : www.nfi.or.th E-mail : cal@nfi.or.th



## Calibration Report

**Certificate No.:** 2202093-001-01  
**Equipment:** Digital Thermometer with RTD (pH Meter)  
**Resolution:** 0.1 °C **Model:** SevenEasy pH  
**Serial No.:** 1230525212 **ID No.:** UAE.WAS.003/2553  
**Manufacturer:** METTLER TOLEDO  
**Date of Calibration:** 16 March 2022

Page 5 of 5

**Calibration point:** 15.0, 25.0 and 35.0 °C

### Calibration result:

- The probe was immersed in liquid bath or dry bath to a minimum depth of 120 mm.
- Description of probe, model : N/A S/N : N/A
- Dimension of probe : Diameter 3.5 mm., Length 135 mm.,
- Sheath material : Stainless Steel

UUC* Reading (°C)	Standard Temperature (°C)	Correction Value (°C)	Uncertainty ± (°C)
15.2	15.001	-0.2	0.099
25.2	25.002	-0.2	0.099
35.2	35.002	-0.2	0.099

**Note** - UUC\* : Unit Under Calibration

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor  $k=2$ , providing a approximately 95 %.

----- End -----

เอกสาร



มูลนิธิส่งเสริมพัฒนาอุตสาหกรรม  
ศูนย์บริการห้องปฏิบัติการอุตสาหกรรมอาหาร

Foundation for Industrial Development National Food Institute  
Food Industrial Laboratory Service Center



## Calibration Certificate

**Certificate No.:** 2203527-001-01  
**Client name:** UNITED ANALYST AND ENGINEERING CONSULTANT CO., LTD.  
**Address:** 3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchack, Prakhnong, Bangkok 10260

Page 1 of 5

**Equipment:** pH Meter  
**Manufacturer:** METTLER TOLEDO  
**Model:** Seven Compact S220  
**Serial No.:** C113432421  
**ID No.:** UAE.WAT.009/2564  
**Order No.:** 2203527  
**Operation No.:** 2203527-001  
**Date of Receipt:** 30 June 2022  
**Date of Calibration:** 5 July 2022

**Calibrated by** Mr.Worapob Sooktong  
Scientist

**Approved by**  
Manager,

**Date of Issue:** 5 July 2022

Responsible for the Technical Management Team

The uncertainties are for a confidence probability of approximately 95%.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65

2008 ซอย 36 ถนนสุขุมวิท แขวงบางค้อ เขตบางพลี กรุงเทพมหานคร เอกสารไม่ควบคุม  
2008 Soi 36, Arun Amarin Road, Bang Yi Khan Subdistrict, Bang Phlat District, Bangkok 10700, Thailand  
Tel: +66(0) 2422 8588 Fax: +66(0) 2422 8545



nfi.or.th

## Calibration Report

**Certificate No.:** 2203527-001-01  
**Equipment:** pH Meter  
**Resolution:** 0.01 pH ; 1 mV  
**Manufacturer:** METTLER TOLEDO  
**Model:** Seven Compact S220  
**Serial No.:** C113432421  
**Type:** Bench top  
**ID No.:** UAE.WAT.009/2564  
**Date of Calibration:** 5 July 2022 Page 2 of 5

**Location:** Calibration Laboratory, National Food Institute  
**Environment Condition:** Ambient Temperature: ( 23.5 ± 1.5 ) °C Relative Humidity: ( 53 ± 5 ) %  
**Condition of Equipment:** Good Condition  
**Condition of this Results of Calibration**

1. Calibration Method In house method : W-CC-002 based on direct measurement by using standard voltage calibrator and certified reference material (CRM)

2. Reference Standards / Certified Reference Material

Instruments	Serial / ID No.	Manufacturer	Certificate No.	Due Date
2.1 DC Voltage Calibrator	2709007	Fluke	22E1959	17 June 2023
2.2 Digital Thermometer	2709007	Fluke	CC-640599-01	30 October 2022
2.3 Thermo-Hygro Meter	NFI.BTH005/18	PONPE	QR22-0351	18 February 2023
Certified Reference Material	Lot No.	Manufacturer	Ref N	Expire Date
2.4 pH buffer 4.008 (Primary pH buffer Solution)	805203	CPAchem	PH216.L5	21 April 2024
2.5 pH buffer 6.865 (Primary pH buffer Solution)	805204	CPAchem	PH217.L5	21 April 2024
2.6 pH buffer 10.01 (Primary pH buffer Solution)	805205	CPAchem	PH220.L5	21 April 2023
2.7 pH buffer 7.00 (Standard pH buffer Solution)	805206	CPAchem	PH107.L5	21 April 2023

3. This certification is traceable to The International System of Unit (SI Unit)

3.1 Instruments No.2.1	through	NSC-TISI-TIS 17025 Laboratory Accreditation of Calibration No.0075
3.2 Instruments No.2.2	through	NSC-TISI-TIS 17025 Laboratory Accreditation of Calibration No.0061
3.3 Instruments No.2.3	through	NSC-TISI-TIS 17025 Laboratory Accreditation of Calibration No.0292
3.4 Certified Reference Material No. 2.4 to 2.6	traceable to	Primary measurement method- Harned cell using calibrated thermometer, barometer, and nanovoltmeter. The Standard Solution preparation and certified by CPAchem Ltd is accredited to ISO 17034 and ISO/IEC 17025
3.5 Certified Reference Material No.2.7	traceable to	BIM RefN HI-27 LotN 04.06.2021; BIM RefN HI-28 LotN 28.05.2021; BIM RefN HI-27 LotN 04.06.2021; BIM RefN HI-28 LotN 28.05.2021, the Standard Solution preparation and certified by CPAchem Ltd is accredited to ISO 17034 and ISO/IEC 17025

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

F-CS-012 Revision: 01 Date: 20-04-65



## Calibration Report

**Certificate No.:** 2203527-001-01  
**Equipment:** pH Meter  
**Resolution:** 0.01 pH ; 1 mV  
**Manufacturer:** METTLER TOLEDO  
**Model:** Seven Compact S220  
**Serial No.:** C113432421  
**Type:** Bench top  
**ID No.:** UAE.WAT.009/2564  
**Date of Calibration:** 5 July 2022 Page 3 of 5

**Calibration Results:**

1. Calibration of pH Meter (Manual Temperature Compensation at 25 °C)

Nominal pH	DC Voltage Standard (mV)	Average Indicator Reading		Uncertainty (±mV)	Coverage Factor (k)
		mV	pH		
0	414.117	414	0.00	0.58	2.00
2	295.811	296	2.00	0.58	2.00
4	177.462	177	4.00	0.58	2.00
6	59.159	59	6.00	0.58	2.00
7	-0.001	0	7.00	0.58	2.00
8	-59.159	-59	8.00	0.58	2.00
10	-177.463	-177	10.00	0.58	2.00
12	-295.812	-296	12.00	0.58	2.00
14	-414.119	-414	14.00	0.58	2.00

2. Calibration of pH Meter with Electrode (Manual Temperature Compensation at 25 °C)

**Equipment:** pH Electrode  
**Type:** Combined Electrode  
**Manufacturer:** METTLER TOLEDO  
**Model:** InLab Expert Pro-ISM  
**Serial No.:** 2210418  
**ID.No.:** N/A

**Performance of Electrode system** (Three-Point Calibration at pH4, pH7 and pH10)

Certified Value @25 °C (pH)	Average Indicator Reading		Relative Slope (%)	Uncertainty (± pH)	Coverage Factor (k)
	pH	mV			
4.008	4.010	182	-	0.0071	2.00
6.865	6.850	14	100.0	0.0075	2.00
10.008	10.010	-169	97.9	0.0093	2.00
6.985	6.990	6	-	0.0087	2.00

F-CS-012 Revision: 01 Date: 20-04-65



## Calibration Report

**Certificate No.:** 2203527-001-01  
**Equipment:** Digital Thermometer with RTD  
 Resolution: 0.1 °C Model: Seven Compact S220  
 Serial No.: C113432421 ID No.: UAE.WAT.009/2564  
 Manufacturer: METTLER TOLEDO  
**Date of Calibration:** 5 July 2022

Page 4 of 5

**Location:** Calibration Laboratory, National Food Institute  
**Environment Condition:** Ambient Temperature 25 °C ± 1 °C  
 Relative Humidity 48 % ± 3 %

### Condition of this results of Calibration:

- Calibration Method : - In house method: W-TE-025 by comparison with standard thermometer.  
 - The Calibration is determined by comparing with a known temperature from a standard resistance thermometer.  
 - The temperature scale in use at this laboratory is the International Temperature scale of 1990 ( ITS-90 ).

### 2. Reference Standard Instrument :

Instrument	Model	Serial No.	Certificate No.	Due Date	Through
HANDHELD THERMOMETER	1521	A85997	TE 650057-01	10-Dec-22	NATIONAL FOOD INSTITUTE
Platinum Resistance Thermometer (PRT)	385	509201			

Support Equipment : - Low Temperature Bath (ISOCAL-6), Model: Europa-6 Plus Basic, S/N: 341592/2

- This certificate is traceable to International System of Units (SI Units).
- This certificate was certified only for the Instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.
- Condition of Calibrated item : Good
- Result of Calibration : ☒ Without adjustment ☐ After adjustment

F-CS-012 Revision: 01 Date: 20-04-65

เอกสารไม่ควบคุม



nfi.or.th

## Calibration Report

**Certificate No.:** 2203527-001-01  
**Equipment:** Digital Thermometer with RTD  
 Resolution: 0.1 °C Model: Seven Compact S220  
 Serial No.: C113432421 ID No.: UAE.WAT.009/2564  
 Manufacturer: METTLER TOLEDO  
**Date of Calibration:** 5 July 2022

Page 5 of 5

**Calibration point:** 15.0, 25.0 and 35.0 °C  
**Calibration result:**

- The probe was immersed in liquid bath or dry bath to a minimum depth of 120 mm.
- Description of probe, model : - S/N : -  
 Dimension of probe : Diameter 9 mm., Length 120 mm.,  
 Sheath material : Stainless Steel

UUC* Reading (°C)	Standard Temperature (°C)	Correction Value (°C)	Uncertainty ± (°C)
15.1	15.038	- 0.1	0.12
25.1	25.038	- 0.1	0.12
35.2	35.024	- 0.2	0.12

### Note

- UUC\* : Unit Under Calibration

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor of approximately 95 %.

----- End -----

F-CS-012 Revision: 01 Date: 20-04-65

เอกสารไม่ควบคุม



nfi.or.th



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES  
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250  
TEL. 0-2717-3000-27 FAX. 0-2719-9484



Cert.No.: 22MM210  
Page.: 1 of 3

## Certificate of Calibration

**Equipment :** Electronic Balance  
**Manufacturer :** Mettler Toledo  
**Model :** XSR205  
**Serial No. :** C009071872  
**ID No. :** UAE.WAO.012/2563  
**Submitted by :** United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phakhanong,  
Bangkok 10260  
**Location :** Balance Room  
**Received order :** 26 April 2022  
**Calibration Date :** 26 April 2022  
**Ambient Temperature :** 15 °C to 40 °C  
**Relative Humidity :** 30 % to 90 %  
**Calibrated by :** Kunchit Promprat  
**Approved by :**   
( ) Pornthippa Tameyakul  
(✓) Malee Butkruea  
( ) Suwit Imjai  
**Issue Date :** 29 April 2022

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



**Equipment :** Electronic Balance  
**Condition As-Received :** Used Item  
**Reference :** 2204-0542OC-1  
**Procedure used :-**

Cert.No.: 22MM210  
Page: 2 of 3

Calibration were conducted using in-house calibration procedure CP-OB01 according to direct measurement method against standard weight.

### Condition of this result of calibration

#### 1. Reference standard instruments:-

Instruments	Model	Serial No.	ID No.	Test report No.	Due date
1) Standard Weight Set (E2)	15884	-	70RC138	MM-0009-21	3 Feb 2023

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This result of calibration was made on requested at the point specified by customer.

4. This certificate is not certified for any commercial transaction.

5. This certification is traceable to the International System of Unit.

**Result of calibration** ( ) Without Adjustment ( \* ) After Adjustment by Internal Calibration

Range capacity :	0 g to 81 g	Resolution	0.00001 g
	81 g to 220 g	Resolution	0.0001 g

#### Before Adjustment :

Applied Weight	Balance Reading	Correction	Measurement Uncertainty	Coverage Factor
( g )	( g )	( g )	( ± mg )	( k )
80	80.00004	-0.00004	0.15	2.00
200	199.9999	+0.0001	0.35	2.00

#### After Adjustment :

#### 1. Determination of the standard deviation of weighing machine ( n = 10 )

Applied Weight	Standard Deviation
( g )	of Reading ( g )
80	0.000008
200	0.00005

เอกสารไม่ควบคุม




Equipment : Electronic Balance  
Condition As-Received : Used Item  
Reference : 2204-0542OC-1  
Result of calibration

## 2. Effect of off center loading

A mass of 100 g was placed to various position on the pan.  
The weighing machine reading error obtained is given in the table

Position 1 (g)	Position 2 (g)	Position 3 (g)	Position 4 (g)	Position 5 (g)
-0.0002	-0.0001	0.0000	-0.0002	-0.0002

Cert.No.: 22MM210  
Page: 3 of 3

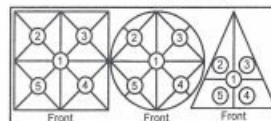

Maximum difference between  
off-center and central loading  
(g)  
0.0002

## 3. Departure from nominal value

Applied Weight (g)	Balance Reading (g)	Correction (g)	Measurement Uncertainty (± mg)	Coverage Factor (k)
Unload	0.00000	0.00000	0.016	2.13
0.05	0.05001	-0.00001	0.016	2.13
0.1	0.10001	-0.00001	0.017	2.11
1	1.00002	-0.00002	0.019	2.05
5	5.00003	-0.00003	0.026	2.00
20	20.00008	-0.00008	0.049	2.00
50	50.00010	-0.00010	0.080	2.00
80	80.00014	-0.00014	0.15	2.00
100	100.0001	-0.0001	0.21	2.00
150	150.0001	-0.0001	0.29	2.00
200	200.0001	-0.0001	0.35	2.00

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor  $k$ , providing a level of confidence of approximately 95 %.

-000-

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES  
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250  
TEL. 0-2717-3000-27 FAX. 0-2719-9484



Cert. No.: 22TM1490  
Page : 1 of 3

# Certificate of Calibration

Equipment : Hot Air Oven  
Manufacturer : Memmert  
Model : UF 55  
Serial No. : B216.1666  
ID No. : UAE.WAO.027/2559  
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udumsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260  
Location : Lab Floor 2  
Received Order : 19 October 2022  
Calibration Date : 19 October 2022  
Ambient Temperature : ( 26 ± 10 ) °C  
Relative Humidity : ( 50 ± 30 ) %

Calibrated by : Preecha Hlahib

Approved by :

Approved Signatory

( ) Ponthippa Tameyakul  
( ) Malee Butkruea  
( ✓ ) Suwit Imjai

Issue Date : 31 October 2022

The Uncertainties are for a confidence probability of approximately 95 %

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม

A 0046800



Equipment : Hot Air Oven  
 Condition As-Received : Used Item  
 Reference : 2210-0575OC-1  
 Procedure Used :-

Cert. No.: 22TM1490  
 Page : 2 of 3

Calibration were conducted using calibration procedure CP-OT02 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector ( RTD ) and Thermocouple Type T.

The temperature scale used was based on ITS-90.

#### Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1 ) Data Acquisition	34970A	MY41021843	22LM4	10 Jan 2023


2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

**Result of Calibration :-** ( \* ) Without Adjustment

**Function of UUC\* :** Temperature Source

**Fresh air setting :** Close

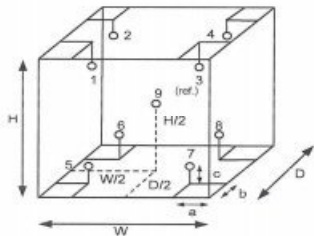

Environment during calibration		
	Beginning	Finished
Temp. ( °C )	29	30
REL.Humid. ( % )	47	40
AC Supply ( Volt )	221	220

Probe Installation Details :		Dimension of Chamber :	
a =	5.0 cm	D =	0.33 m
b =	5.0 cm	W =	0.40 m
c =	5.0 cm	H =	0.40 m
		Capacity =	0.053 m <sup>3</sup>

Ref. Std. ID No.: @ Calibration Point		
Position :	( 104 ) °C	( 140,180 ) °C
1	18-04RTD-01	21-04TC-01
2	18-04RTD-02	21-04TC-02
3	18-04RTD-03	21-04TC-03
4	18-04RTD-04	21-04TC-04
5	18-04RTD-05	21-04TC-05
6	18-04RTD-06	21-04TC-06
7	18-04RTD-07	21-04TC-07
8	18-04RTD-08	21-04TC-08
9 (ref.)	18-04RTD-09	21-04TC-09

เอกสารไม่คว

a 1133252



Equipment : Hot Air Oven  
 Condition As-Received : Used Item  
 Reference : 2210-0575OC-1  
 Result of Calibration :- ( \* ) Without Adjustment

Cert. No.: 22TM1490  
 Page : 3 of 3

**Function of UUC\* :** Temperature Source

**Fresh air setting :** Close

Calibration Point ( °C )	UUC* Setting ( °C )	UUC* Reading ( °C )	Temperature stability ( ± °C )	Temperature uniformity ( °C )	Overall Variation ( °C )	Uncertainty ( ± °C )	Coverage Factor k
104.0	104.0	104.0	0.061	1.3	1.7	0.42	2
140.0	140.0	140.0	0.14	2.3	2.4	1.1	2
180.0	180.0	180.0	0.21	3.5	3.6	1.3	2

Calibration Point ( °C )	Measured Temperature ( °C )								
	Position								
	1	2	3	4	5	6	7	8	9 (ref.)
104.0	103.076	103.876	103.777	104.124	104.667	104.426	104.012	103.928	104.370
140.0	138.199	139.189	138.808	139.550	140.266	139.622	139.293	139.385	140.369
180.0	177.930	179.267	178.643	179.753	181.011	180.093	179.496	179.743	181.278

**Average\* :** The average of 30 values in each position.

**Temperature stability :** One-half of the greatest maximum difference of measured temperature at any one sensor

**Temperature uniformity :** The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

**Overall Variation :** The Difference of the maximum and minimum measured temperatures throughout observation

**UUC\* :** Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่คว

a 1133251

## Calibration Certificate

**Certificate No.:** 2202934-001-01  
**Client name:** UNITED ANALYST AND ENGINEERING CONSULTANT CO., LTD.  
**Address:** 3 Soi Udomsuk 41, Sukhumvit Road,  
 Bangchack, Prakhong, Bangkok 10260

Page 1 of 4

**Equipment:** Electronic Balance

**Manufacturer:** METTLER TOLEDO

**Model:** XSR204

**Serial No.:** C117635043

**ID No.:** UAE.WAS.012/2564


**Order No.:** 2202934

**Operation No.:** 2202934-001

**Date of Receipt:** 13 May 2022

**Date of Calibration:** 13 May 2022

**Calibrated by** Mr.Manas Somsak  
 Specialist

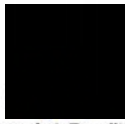**Approved by**   
 ( Mr.Pheraphat Tuanjit )  
 Manager, Division of Calibration Laboratory  
 Responsible for the Technical Management Team

**Date of Issue:** 25 May 2022

The uncertainties are for a confidence probability of approximately 95%

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65

## Calibration Report

**Certificate No.:** 2202934-001-01  
**Equipment:** Electronic Balance **Manufacturer:** METTLER TOLEDO  
**Model:** XSR204 **Resolution:** 0.0001 g  
**Serial No.:** C117635043 **ID No.:** UAE.WAS.012/2564  
**Capacity:** 220 g

**Date of Calibration:** 13 May 2022

Page 2 of 4

**Environment Condition:** Ambient Temperature 22.3 ± 0.1 °C Relative Humidity: 47 ± 3 %

**Place of Calibration:** Balance room (Water Analysis Unit), UNITED ANALYST AND ENGINEERING CONSULTANT CO., LTD.

**Condition of Equipment:** Good Condition

**Condition of This Results of Calibration:**

1. Calibration Method: NFI Method W-MA-001 In-House Method based on UKAS Lab 14 : 2019
2. Reference Standards:

Reference Standard	Model	Serial No	Calibrated By	Certificate No.	Due Date
Standard Weight Class E2	1mg to 200g	B505567572	TCS	M22041375	23 April 2023

Instrument	Model	Serial No	Calibrated By	Certificate No.	Due Date
Thermo-Hygro Meter	PONPE 490	NFL.BTH 010/18	Quality Reborn	QR22-0350	18 February 2023

3. This certification is traceable to SI UNIT

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

**Calibration Results:**


**1. Repeatability of Reading:**

Nominal Value ( g )	Standard Deviation of Reading ( g )
100	0.000033
200	0.000032

**2. Off-Center Error:**

A mass of 50 g was placed and moved to various position on pan.

The balance reading obtained is given in the table.

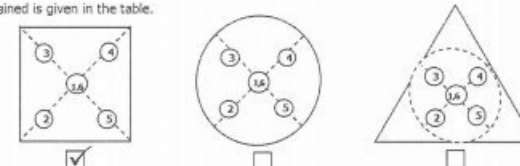

1	2	3	4	5	6	(Maximum Difference)
( g )	( g )	( g )	( g )	( g )	( g )	( g )
50.000	50.000	50.000	50.000	50.000	50.000	0.000

F-CS-012 Revision: 01 Date: 20-04-65

## Calibration Report

**Certificate No.:** 2202934-001-01

**Equipment:**

Electronic Balance

**Manufacturer:** METTLER TOLEDO

**Model:** XSR204

**Resolution:** 0.0001 g

**Serial No.:** C117635043

**ID No.:** UAE.WAS.012/2564

**Capacity:** 220 g

**Date of Calibration:** 13 May 2022

Page 3 of 4

**Calibration Results:** (Continued)

**Calibration Range:** 0 - 200 g

**Calibration Adjustment:** Internal Calibration

**3. Departure from Nominal Value:**

Nominal Value ( g )	Standard Value ( g )	Average Reading ( g )	Correction ( g )	Uncertainty ( ± g )	Coverage Factor k
Unload	0.00000	0.0000	0.0000	0.000085	2.00
0.01	0.01000	0.0100	0.0000	0.000085	2.00
0.02	0.02000	0.0200	0.0000	0.000085	2.00
0.05	0.05000	0.0500	0.0000	0.000085	2.00
0.1	0.10001	0.1000	0.0000	0.000085	2.00
0.2	0.20001	0.2000	0.0000	0.000085	2.00
0.5	0.50002	0.5000	0.0000	0.000085	2.00
1	1.00001	1.0000	0.0000	0.000086	2.00
2	2.00003	2.0000	0.0000	0.000086	2.00
3	3.00004	3.0000	0.0000	0.000087	2.00
5	5.00002	5.0000	0.0000	0.000087	2.00
10	10.00001	10.0000	0.0000	0.000088	2.00
20	20.00004	20.0000	0.0000	0.000092	2.00
30	30.00005	30.0001	-0.0001	0.00010	2.00
40	40.00008	40.0001	0.0000	0.00011	2.00
45	45.00010	45.0001	0.0000	0.00013	2.00

F-CS-012 Revision: 01 Date: 20-04-65

เอกสารไม่ควบคุม



nfi.orth

## Calibration Report

**Certificate No.:** 2202934-001-01

**Equipment:**

Electronic Balance

**Manufacturer:** METTLER TOLEDO

**Model:** XSR204

**Resolution:** 0.0001 g

**Serial No.:** C117635043

**ID No.:** UAE.WAS.012/2564

**Capacity:** 220 g

**Date of Calibration:** 13 May 2022

Page 4 of 4

**Calibration Results:** (Continued)

**Calibration Range:** 0 - 200 g

**Calibration Adjustment:** Internal Calibration

**3. Departure from Nominal Value:**

Nominal Value ( g )	Standard Value ( g )	Average Reading ( g )	Correction ( g )	Uncertainty ( ± g )	Coverage Factor k
50	50.00004	50.0001	-0.0001	0.00011	2.00
55	55.00006	55.0001	0.0000	0.00012	2.00
60	60.00005	60.0001	-0.0001	0.00012	2.00
65	65.00007	65.0002	-0.0001	0.00013	2.00
70	70.00008	70.0002	-0.0001	0.00013	2.00
75	75.00010	75.0002	-0.0001	0.00013	2.00
80	80.00009	80.0002	-0.0001	0.00014	2.00
85	85.00011	85.0002	-0.0001	0.00014	2.00
90	90.00012	90.0002	-0.0001	0.00015	2.00
100	100.00008	100.0003	-0.0002	0.00016	2.00
120	120.00011	120.0003	-0.0002	0.00018	2.00
150	150.00012	150.0004	-0.0003	0.00021	2.00
200	200.00015	200.0004	-0.0003	0.00028	2.00

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor providing a level of confidence of approximately 95 %.

----- End -----

F-CS-012 Revision: 01 Date: 20-04-65

เอกสารไม่ควบคุม



nfi.orth



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES  
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250  
TEL. 0-2717-3000-27 FAX. 0-2719-9484



## Certificate of Calibration

Cert. No.: 22TM90  
Page.: 1 of 3

**Equipment :** BOD Incubator  
**Manufacturer :** Arco  
**Model :** UC4-1320  
**Serial No. :** 13URC4S013201  
**ID No. :** UAE.WAO.015/2561  
**Submitted by :** United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260  
**Location :** Lab Floor 2  
**Received Order :** 17 February 2022  
**Calibration Date :** 17 February 2022  
**Ambient Temperature :**  $(26 \pm 10) ^\circ\text{C}$   
**Relative Humidity :**  $(50 \pm 30) \%$

**Calibrated by :**

**Approved by :**

( ) Pornthippa Tameyakul  
( ) Malee Butkruea  
( ) Suwit Imjai

**Issue Date :** 22 February 2022

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม  
A 0038099



**Equipment :** BOD Incubator  
**Condition As-Received :** Used Item  
**Reference :** 2202-0446OC-1  
**Procedure Used :-**

Cert. No.: 22TM90  
Page.: 2 of 3

Calibration were conducted using calibration procedure CP-OT02 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector ( RTD ).

The temperature scale used was based on ITS-90.

### Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1 ) Data Acquisition	34970A	MY44035217	21LM30	23 Dec 2022


2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

**Result of Calibration :-** ( \* ) Without Adjustment

**Function of UUC\* :** Temperature Source

**Fresh air setting :** Not Available

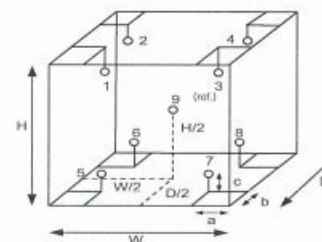

**Probe Installation Details :**

a = 10 cm  
b = 10 cm  
c = 10 cm

**Dimension of Chamber :**

D = 0.62 m  
W = 1.2 m  
H = 1.2 m  
Capacity = 0.89 m<sup>3</sup>

Environment during calibration		
	Beginning	Finished
Temp. ( °C )	28	28
REL.Humid. ( % )	68	75
AC Supply ( Volt )	226	226

Position :	Ref. Std. ID No.:
1	18-10RTD-01
2	18-10RTD-02
3	18-10RTD-03
4	18-10RTD-04
5	18-10RTD-05
6	22-10RTD-10
7	18-10RTD-07
8	18-10RTD-08
9 (ref.)	18-10RTD-09

เอกสารไม่ควบคุม  
a 1096042



Equipment : BOD Incubator  
Condition As-Received : Used Item  
Reference : 2202-0446OC-1  
Result of Calibration :- ( \* ) Without Adjustment  
Function of UUC\* : Temperature Source  
Fresh air setting : Not Available

Cert. No.: 22TM90  
Page.: 3 of 3

Calibration Point ( °C )	UUC* Setting ( °C )	UUC* Reading ( °C )	Temperature stability ( ± °C )	Temperature uniformity ( °C )	Overall Variation ( °C )	Uncertainty ( ± °C )	Coverage Factor k
20.0	19.5	19.4	0.30	0.58	1.0	0.55	2

Calibration Point ( °C )	Measured Temperature ( °C )								
	Position								
	1	2	3	4	5	6	7	8	9 (ref.)
20.0	20.154	20.013	20.356	19.939	19.834	19.761	19.817	19.824	19.922

Average\* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC\* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่ควบคุม  
a 1096041



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES  
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250  
TEL. 0-2717-3000-27 FAX. 0-2719-9484



Cert. No.: 22TM305  
Page.: 1 of 3

## Certificate of Calibration

Equipment : BOD Incubator  
Manufacturer : ARCO  
Model : UR-1320  
Serial No. : -  
ID No. : UAE.WAO.018/2551  
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260  
Location : Lab Floor 2  
Received Order : 7 April 2022  
Calibration Date : 7 April 2022  
Ambient Temperature : ( 26 ± 10 ) °C  
Relative Humidity : ( 50 ± 30 ) %

Calibrated by : Man Pattanapongpaiboon

Approved by :

( ) Pornthippa Tameyakul  
( ) Malee Butkruea  
( ) Suwit Imjai

Issue Date : 18 April 2022

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม  
A 0040246



Equipment : BOD Incubator  
Condition As-Received : Used Item  
Reference : 2204-00150C-2

Cert. No.: 22TM305  
Page.: 2 of 3

**Procedure Used :-**

Calibration were conducted using calibration procedure CP-OT02 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector ( RTD ).

The temperature scale used was based on ITS-90.

**Condition of this result of calibration**

1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1 ) Data Acquisition	34970A	MY41021843	22LM4	10 Jan 2023


2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

**Result of Calibration :-** ( \* ) Without Adjustment

**Function of UUC\* :** Temperature Source

**Fresh air setting :** Not Available

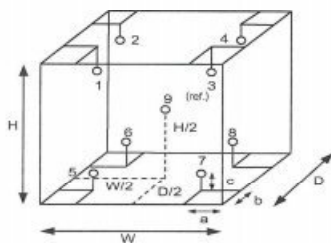

**Probe Installation Details :**

a = 10 cm  
b = 10 cm  
c = 10 cm

**Dimension of Chamber :**

D = 0.62 m  
W = 1.2 m  
H = 1.2 m  
Capacity = 0.89 m<sup>3</sup>

Environment during calibration		
	Beginning	Finished
Temp. ( °C )	27	27
REL.Humid. ( % )	56	59
AC Supply ( Volt )	222	221

Position :	Ref. Std. ID No.:
1	18-04RTD-01
2	18-04RTD-02
3	18-04RTD-03
4	18-04RTD-04
5	18-04RTD-05
6	18-04RTD-06
7	18-04RTD-07
8	18-04RTD-08
9 (ref.)	18-04RTD-09

เอกสารไม่ครบชุด

a 1104314



Equipment : BOD Incubator  
Condition As-Received : Used Item  
Reference : 2204-00150C-2

Cert. No.: 22TM305  
Page.: 3 of 3

**Result of Calibration :-** ( \* ) Without Adjustment

**Function of UUC\* :** Temperature Source

**Fresh air setting :** Not Available

Calibration Point ( °C )	UUC* Setting ( °C )	UUC* Reading ( °C )	Temperature stability ( ± °C )	Temperature uniformity ( °C )	Overall Variation ( °C )	Uncertainty ( ± °C )	Coverage Factor k
20.0	20.0	20.0	0.50	0.44	1.1	0.64	2

Calibration Point ( °C )	Measured Temperature ( °C )								
	Position								
	1	2	3	4	5	6	7	8	9 (ref.)
20.0	20.080	20.056	19.866	19.826	19.655	19.656	19.819	19.979	19.899

**Average\* :** The average of 30 values in each position.

**Temperature stability :** One-half of the greatest maximum difference of measured temperature at any one sensor

**Temperature uniformity :** The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

**Overall Variation :** The Difference of the maximum and minimum measured temperatures throughout observation.

**UUC\* :** Unit Under Calibration

**Note :** The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่ครบชุด

a 1104313

## CERTIFICATE OF CALIBRATION

**Equipment :** COD Test Tube Heater  
**Meter Model :** HI839800-02 **Serial No. :** H0185001  
**Manufacturer :** Hanna Instruments  
**Made in :** Romania  
**Condition As-Received :** Used Product  
**Reference :** RE220234  
**Customer name :** United Analyst and Engineering Consultant Co., Ltd.  
3 Soi Udomsuk 41, Sukhumvit Rd., Bangchak,  
Phrakhanong, Bangkok 10260  
**Received date :** 21 February 2022  
**Calibrate date :** 1 March 2022  
**Issue date :** 2 March 2022  
**Ambient Temperature :**  $(25 \pm 2) ^\circ\text{C}$   
**Relative Humidity :**  $(50 \pm 15) \% \text{ RH}$   
**Calibrated Location :** Hanna Instruments (Thailand) Ltd.

Calibrated by :

Mr. Pichit Petthong  
Calibration Engineer

Approved by :

Mr. Anan Suwanchaisakul  
Authorized Signatory

This certificate was certified only for the instrument we calibrated.

This result of calibration was found accurate on date and place of calibration only.

\*\* This certificate may not be reproduced other than in full, except with the prior written \*\*  
approval of the head of Hanna Instrument (Thailand)

เอกสารไม่ควบคุม

## Condition of this result of calibration

## Reference Standard Instruments :

Instruments	Model	Serial No.	Certificate No.	Traceable
Thermometer With Sensor	HI935005	03250060101	21T167	Technology Promotion Association ( Thailand-Japan )

## Reference / Procedure :

This equipment was calibration by comparison to the reference standard (Standard platinum resistance thermometer) whose accuracy is traceable to the national standard. The calibration was performed by generating the specified working point of temperature then recorded the temperature reading values against the reference standard according to Hanna Calibration Laboratory work Instruction No. 141.

This temperature scale used was based on ITS-90

All data shown below were as-received values without adjustment.

## SITE CALIBRATION

1	2	3	4	5
6	7	8	9	10
11	12	13	14	15
16	17	18	19	20
21	22	23	24	25

เอกสารไม่ควบคุม

**Result of Calibration :**

Calibration Point	Unit Under Calibration Setting	Unit Under Calibration Reading	Temperature Stability	Uncertainty of Measurement
150.0 (°C)	- (°C)	150.6 (°C)	1.3 (°C)	± 0.39 (°C)

Calibration Point (°C)	Average Standard Reading (°C)				
	Position				
150.0	1	2	3	4	5
	150.2	150.4	150.4	150.3	150.2
	6	7	8	9	10
	150.4	150.9	151.1	151.1	150.6
	11	12	13	14	15
	150.4	151.0	151.5	151.3	150.5
	16	17	18	19	20
	150.3	150.8	151.2	151.2	150.5
	21	22	23	24	25
	150.2	150.3	150.5	150.4	150.3

The reported uncertainty of measurement was based on standard uncertainty multiplied by a coverage factor k=2, providing a level of confidence of approximately 95%

\*\* End of certificate \*\*

เอกสารไม่ควบคุม

**Verification Certificate**

**Certificate No.:** 2202361-001-01  
**Client name:** UNITED ANALYST AND ENGINEERING CONSULTANT CO., LTD.  
**Address:** 3 Soi Udomsuk 41, Sukhumvit Road, Bangchack, Prakhong, Bangkok 10260

Page 1 of 4

**Equipment:** HEATING BLOCK DIGESTION

**Manufacturer:** FOSS

**Model:** 2520

**Serial No.:** 91794469

**ID No.:** UAE.WAS.011/2560

**Order No.:** 2202361

**Operation No.:** 2202361-001

**Date of Receipt:** 4 April 2022

**Date of Calibration:** 4-6 April 2022

**Calibrated by** Mr.Nuttapol Niyomchat  
Specialist

**Approved by**

Manager, Division of Calibration Laboratory  
Responsible for the Technical Management Team

**Date of Issue:** 11 April 2022

The uncertainties are for a confidence probability of approximately 95 %.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

เอกสารไม่ควบคุม

## Verification Report

**Certificate No.:** 2202361-001-01  
**Equipment:** HEATING BLOCK DIGESTION  
 Model: 2520 Serial No.: 91794469  
 Resolution: 1 °C ID No.: UAE.WAS.011/2560  
 Manufacturer: FOSS  
**Date of Calibration:** 4-6 April 2022

Page 2 of 4

**Location:** Laboratory Room, NATIONAL FOOD INSTITUTE  
**Environment Condition:** Ambient Temperature ( 25 ± 3 ) °C  
 Relative Humidity ( 55 ± 15 ) %  
 Line Voltage ( 220 ± 10 ) Volt

## Condition of this results of Calibration:

- This instrument was calibrated by insert standard thermocouples type R into its heating block digestion and compared to temperature obtained from reference standards thermometer at calibrated point.  
 - The temperature scale used was based on ITS - 90 .  
 - All data show below were final values and the initial data may be obtained upon request.

## 2. Reference Standard Instrument :

Instrument	Model	Serial No.	Certificate No.	Due Date	Through
Digital Thermometer with Thermocouple	34970A/34901A	MY44045576/MY41194453	TC21/0041	24-Apr-2022	N.M. Technical Center Laboratory
	Type R	TC#101-103 / CH#101-103			

- This certificate is traceable to international system of units (SI Units).
- This certificate was certified only for the instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.
- Condition of Calibrated item : Good

## UUC\* Description

Time of Record - Hour 30 Minute At 380 °C

7. Result of Calibration : ☒ Without adjustment ☐ After adjustment

## Verification Report

**Certificate No.:** 2202361-001-01  
**Equipment:** HEATING BLOCK DIGESTION  
 Model: 2520 Serial No.: 91794469  
 Resolution: 1 °C ID No.: UAE.WAS.011/2560  
 Manufacturer: FOSS  
**Date of Calibration:** 4-6 April 2022

Page 3 of 4

**Calibration point:** 380 °C**Calibration result:**

## Reporting of Temperature

Block No.	UUC* Setting (°C)	UUC* Reading (°C)	Stability (±°C)	Standard Thermometer (°C)	Uncertainty (±°C)
1	380	380	0.13	376.48	1.5
2	380	380	0.12	376.58	1.5
3	380	380	0.12	376.51	1.5
4	380	380	0.14	376.70	1.6
5	380	380	0.18	376.81	1.6
6	380	380	0.12	377.23	1.6
7	380	380	0.12	377.37	1.5
8	380	380	0.13	376.68	1.5
9	380	380	0.14	376.72	1.5
10	380	380	0.18	378.97	1.6
11	380	380	0.25	378.79	1.6
12	380	380	0.11	377.14	1.6
13	380	380	0.19	379.65	1.6
14	380	380	0.16	379.61	1.6
15	380	380	0.16	378.66	1.6
16	380	380	0.15	379.18	1.6
17	380	380	0.23	377.39	1.6
18	380	380	0.11	377.71	1.6
19	380	380	0.22	376.64	1.6
20	380	380	0.16	376.56	1.6

## Note:

- UUC\* = Unit Under Calibration
- Immersion depth of standard thermometer in tube level high of sand is equal heater plate of UUC.
- Stability = One-half of the greatest maximum difference of measured temperatures at one sensors, for at least half an hour after reaching steady state.



## Customer Service Report

FOSS South East Asia  
3388 Sirinrat Building, 25th – 26th Floor, Unit No. 3388/90,  
Rama IV Road, Klongton , Klongtoey, Bangkok, Thailand 10110

Report No:	6534
------------	------

Date:	25/1/2022		Phone	
Customer:	United Analyst and Engineering		Address:	
Instrument:	KT6100		Serial: 91829052	

Hours Start Finish	Travel To Customer		Labour		Travel From Customer	
	2.00	0.5	2-12 am	6 hr	4-6 c	1 hr
	2.30		1-4 pm		5.30	

Job Type					
Application		Special		Standard	
Normal	1	Courtesy Visit	1	Installation	1
Distributor	1	PMA Onboarding	1	Quote	1
Internal	1	Warranty	1	Repair	1
Digital Service	1	Sales Support	1	Remote	1
				Training	1
				In House	1
				PM	1
				Other	1

<b>PO/Quote Number:</b>	0000000000
-------------------------	------------

PMA Type	Applicable	Contract No.	Application
----------	------------	--------------	-------------

Details of Work / Test		Condition / Status
		Pack
Run the software	- Software - program Editor - mp setting - Manual fun - user maintenance - Run Blank - Run Recovery	
Instrument Ready for Use	(OK)	Not OK

[illegible]

I confirm this report is accurate and complete

Signature: [Redacted]

Name: [Redacted]

Would you be willing to participate in a brief survey in order to tell us how we performed? [Redacted]

## เอกสารไม่ควบคุม

**FOSS**

## Installation Qualification

## Kjeltec™ 8100 Distillation Unit

**This IQ applies to Kjeltex™ 8100 Distillation Unit manufactured by FOSS Analytical. The installation is performed by FOSS trained service personnel.**

## 1 Intended Use

Kjeltec 8100 is intended for laboratory use analyzing parameters as specified in FOSS Analytical AB's Application Notes.

## 2 Purpose

This installation Qualification is designed to assure that:

- The Kjeltec instrument is received complete, with all required parts in good condition.
- The location of the instrument is environmentally and ergonomically suitable
- The instrument is assembled and configured correctly
- Suitable electricity and water are supplied to the instrument. see table 2 for requirements.

### 3 Identification

Description	Serial Number
Kjeltec 8100 Distillation Unit	91889052

Dedicated Analytical Solutions

**FOSS Analytical A/S**  
69 Slangerupgade  
DK-3400 Hillerød  
Denmark  
Tel +45 7010 3370  
Fax +45 7010 3371  
E-mail support@foss.dk  
Web www.foss.dk

Customer Support, 6003 7242 / Rev. 1

**FOSS Analytical AB**  
Box 70  
SE-263 21 Högåns  
Sweden

เอกสารไม่ควบคุม

1(6)

## 4 Control of Received Equipment

### 4.1 Verify that the correct instrument type and accessory kit items are received and in proper condition

The packing list (shipped with the instrument) specifies all the items. The installer will verify that all items are received as shipped on the packing list. For each item listed, verify that the acceptance criteria are met. If so, write "Y" in the right column of the table immediately following.

Packing List Item	Acceptance Criteria	Pass(Y/N)
Kjeltec 8100 Distillation Unit	No visible damage, received in undamaged FOSS Analytical's standard shipping container	Y
Accessory kit, according to packing list	Included. No visible damage, received in undamaged FOSS Analytical's standard shipping container	Y
Handling device for digestion tube	Included. No visible damage.	Y
Tanks with level sensors for Waste, Alkali and Water	Included. No visible damage.	Y
Receiver flask	Included. No visible damage.	Y
One digestion tube 250ml One digestion tube 100 ml	Included. No visible damage.	Y
Tube adapter	Included. No visible damage.	Y
User manual	Kjeltec 8100 Distillation Unit	Y
Owners guide	Kjeltec 8100 Distillation Unit	Y
Quick guide	Kjeltec 8100 Distillation Unit	Y
Spare parts manual	Kjeltec 8100 Distillation Unit	Y
Application notes	AN 300 included AN 303 included	Y

## 5 Installation


### 5.1 The equipment must be installed in a suitable location with power, water and draining available

Verify that the instrument installation site meets the acceptance criteria given in the table below. If so, write "Y" in the right column of the table immediately following.

Location Requirements	Acceptance Criteria	Pass (Y/N)
Adequate space for instrument	Dimensions 48x58x69 cm	Y
AC supply available for instrument	200-240 V 50/60Hz	Y
Current	10 A	Y
Cold water supply available	2 L/min at 30°C	Y
Drain	For cooling water and waste (depending on local waste disposal legislation)	Y
Ambient temperature	Max. 40°C	Y
Ambient humidity	Max. 80% relative	Y
Internal fuses	T10A AH	Y

## 5.2 The instrument must be assembled correctly

Verify that all tubes are correct connected. If so, write "Y" in the right column of the table immediately following.

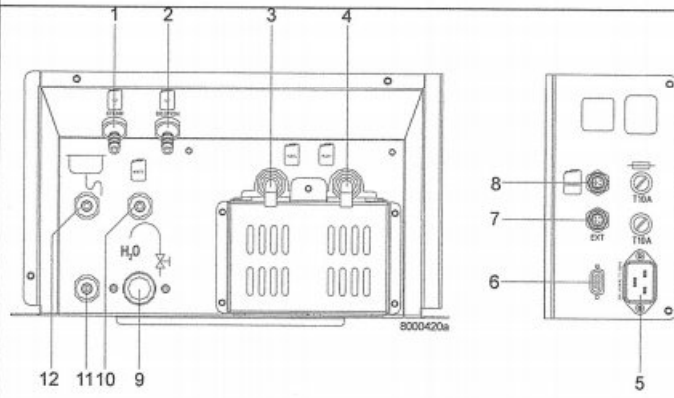Instrument Tubing Connections	Acceptance Criteria	Pass (Y/N)
 <p>1. Deionised water in (steam generator)  2. Deionised water in (dilution water)  3. *) Receiver solution in  4. Alkali in  5. Power  6. Not used  7. External titration module  8. Level sensors  9. Cooling water in (tap water)  10. Waste water out (tube drain vessel)  11. Drain  12. Cooling water out (tap water)  *) Only on Kjeltac 8200</p>	Visual verification by installer	Y

## 5.3 The instrument should be assembled and powered up

Connect the distilling unit to the power supply. Perform the start up procedure and check that the expected response is obtained. If so, write "Y" in the right column of the table immediately following.

Action	Expected Response	Pass (Y/N)
Switch on the power	The instruments start up and the self test will run. The sample counter shows the number of analysed samples since first power and the Software Version shows the version of the instruments software.	Y
	After start-up, Program 1 is loaded and the Analyse menu is displayed.	Y
Turn on the cold water tap	No visible reaction	Y
Press the "Manual" view	The Manual menu is opened	Y
Open the door with the handle, place the test tube and receiver flask in position. Close the door.		Y
Select <b>Dilution</b> and press <b>Start</b>	Water is added to the tube	Y
Select <b>Alkali</b> and press <b>Start</b>	Alkali is added to the tube	Y
Select <b>Steam</b> and press start	After heating up, steam is entering the tube	Y
Select <b>Drain</b> and press <b>Start</b>	The tube is drained	Y

## 6 Summary of Deviations/Comments

Deviations from above requirements are specified below and any corrective actions are noted.

Deviation	Action	Comment

## 7 IQ Documentation

Upon successful completion and recording of all instructions above, sign and date this sheet below.  
If required by customer, leave one signed copy with instrument.

If customer's internal procedures require further reporting or witnessing of results, execute those procedures as required.

Installed By: Pannipa Onnom

Company: Foss SEA

Customer Name: United Analyst and Engineering

Company: United Analyst and Engineering

Date completed: July 25, 2022

## Kjeltec™ 8100 Distillation Unit

This OQ applies to Kjeltec 8100 Distillation Unit manufactured by FOSS Analytical. The operation qualification is performed by FOSS trained service personnel.

### 1 Intended Use

Kjeltec 8100 is intended for laboratory use analyzing parameters as specified in FOSS Analytical Application Notes.

### 2 Purpose

This procedure is designed to test the function of the instrument according to factory test specifications:

- Alkali volume
- Distillation Accuracy
- Distillation Repeatability

### 3 Identification

Description	Serial Number
Kjeltec 8100 Distillation Unit, 200-240 V 50/60 Hz	91899052

#### Dedicated Analytical Solutions

FOSS Analytical A/S  
69 Slangerupgade  
DK-3400 Hillerød  
Denmark

Tel +45 7010 3370  
Fax +45 7010 3371  
E-mail support@foss.dk  
Web www.foss.dk

FOSS Analytical AB  
Box 70  
SE-263 21 Höganäs  
Sweden

Tel +46 42 361500  
Fax +46 42 340349  
E-mail support@foss.dk  
Web www.foss.dk

## 4 Performance

### 4.1 Verify the dispensed volumes of reagents

**Note!** To verify the dispensed volumes of reagents a triple test should be done to be statistic correct. Then calculate a mean value.

1. Choose "Manual" in the menu. (When starting up the instrument Program 1 is loaded)
2. Open the safety door by pressing **Open** and place a tube in the instrument. Close the safety door.

#### Water

1. Press **Dilution** and then press **Start**. 80 ml of water will be filled into the tube.
2. Measure the collected water in a graduated measuring glass and note the result in table 1 below.
3. Check acceptance criteria in the table and make the judgment if passed or not.

**Note!** If the water volume needs to be calibrated, go to 4.8.5 Dilution Pump Calibration in the User Manual.

#### Alkali

1. Press **Alkali** and then press **Start**. 50 ml of alkali will be filled into the tube.
2. Measure the collected alkali in a graduated measuring glass and note the result in table 1 below.
3. Check acceptance criteria in the table and make the judgment if passed or not.

Table 1 Volume control

Test	Result	Expected result	Passed (Y/N)
Water volume	<u>83</u> ml <u>83</u> ml <u>82</u> ml Mean <u>82.67</u> ml	76- 84 ml	Y
Alkali volume	<u>52</u> ml <u>52</u> ml <u>53</u> ml Mean <u>52.33</u> ml	47- 54 ml	Y

### 4.2 Verify the distillation procedure, accuracy and precision

The distillation principle is to convert ammonium ( $\text{NH}_4^+$ ) into ammonia ( $\text{NH}_3$ ) by using an alkali ( $\text{NaOH}$ ) and thereafter steam distil it into a receiver flask containing boric acid and titrate with standard acid solution using colorimetric end-point detection. Ammonium sulphate, a substance with known ammonia content, can be used to check the accuracy of the distillation. The recovery is calculated from obtained result.

The way to perform this test will be described in the following.

#### Chemical Check

Use ammonium sulphate ( $(\text{NH}_4)_2\text{SO}_4$ , purity > 99.5 % \*)

Mol. weight = 132.14 g/mol, Nitrogen content in ammonium sulphate (99.5 %) = 21.09% \*)

#### Analysis conditions according to AN 300

Water	80 ml
Alkali	50 ml NaOH (40%w/w)
Receiver solution	30 ml boric acid (4%)
Distillation time	5 minutes
SAFE	5 seconds
Titrant	0.2N HCl

#### For reagent preparation see Appendix A

1. Start the instrument and run two blanks without chemicals according to above analysis conditions, distil into a receiver flask containing boric acid. Titrate with a standard acid solution using colorimetric end-point detection. If the blanks are less than 0,2 ml continue with the recovery tests:
2. Weigh 0.15 g ammonium sulphate into a tube. Prepare 6 samples (tubes).
3. Run the six samples according to above analysis conditions. Titrate with a standard acid solution using colorimetric end-point detection.
4. Calculate the recovery according to below equations. Expected results of recovery should be 100%±1%.

Recovery test	Result	Expected result	Passed (Y/N)
Blank value (water blank)	1. <u>0.08</u> ml 2. <u>0.14</u> ml	0.05-0.20 ml	Y
Recovery	1. <u>100.30</u> % 2. <u>100.30</u> % 3. <u>100.63</u> % 4. <u>99.01</u> % 5. <u>99.97</u> % 6. <u>100.01</u> %		
Accuracy	Mean Value: <u>100.03</u>	99-101%	Y
Precision	SD: <u>0.557</u>	SD <1%	Y

\*) **Note!** Please also note that the below calculations must be adjusted if other purity levels of ammonium salts are used. A certificate for the chemical supplier should be available

Purity	Nitrogen content
99,5%	21.09%
99,6%	21.12%
99,7%	21.14%
99,8%	21.16%
99,9%	21.18%

$$\% \text{ Nitrogen} = \frac{(ml_{\text{sample}} - ml_{\text{blank}}) \times N \times 14,007 \times 100}{mg_{\text{sample}}} \quad \begin{matrix} \nearrow 0.1095 \\ 21.72 \end{matrix}$$

$N$  = Normality of titrant to 4 places of decimal.

$$\% \text{ Recovery} = \frac{\% \text{ Nitrogen}}{21.09} \times 100$$

mg sample

① 0.1592 23.56

②

③

④

⑤

⑥

## 5 Summary of Deviations/Comments

Deviations from above requirements are specified below and any corrective actions are noted.

Deviation	Action	Comment

## 6 OQ Documentation

Upon successful completion of tests above, sign and date this sheet below. If required by customer, leave one signed copy with instrument.

If customer's internal procedures require further reporting or witnessing of results, execute those procedures as required.

Performed By: \_\_\_\_\_

Company: \_\_\_\_\_

Customer Name: \_\_\_\_\_

Company: \_\_\_\_\_

Date completed: \_\_\_\_\_

## 7 Appendix A

### 7.1 Preparation of Reagents

#### 7.1.1 Alkali

To convert ammonium into ammonia an excess of sodium hydroxide is necessary.

Use 400 g NaOH per litre of solution. Commercially available in concentrations up to 50 %. Do not use concentrations above 40 % as this will lead to crystal formation impairing the function of the pumps. If you can only buy concentrations > 40 %, dilute it before use.

#### 7.1.2 Titrant acid, determination of concentration

To be able to achieve accurate nitrogen / protein results, one must be quite sure that the HCl (hydrochloric acid) concentration is what it is supposed to be. A titration against a predetermined solution of sodium carbonate as described below is thus necessary. Incorrect HCl concentration can otherwise cause substantial errors.

- **Standard substance**

Weigh approx. 10 g of anhydrous sodium carbonate ( $\text{Na}_2\text{CO}_3$ ). Use a mortar to make a fine powder. Dry it for 1 h at 265 °C or 2 h at 200 °C. After cooling in a desiccator, transfer the sodium carbonate to a beaker with a tight lid. Store it in a desiccator.

- **Indicator solutions**

Dissolve 0.1 g methyl red in 100 ml methanol. Dissolve 0.1 g bromocresol green in 100 ml methanol.

- **Procedure**

Weigh approx. 0.4 g of the standard substance, using an analytical balance, note the weight ( $W_1$ ). Transfer the sodium carbonate to a receiver flask and add 40 ml of  $\text{H}_2\text{O}$  (distilled or deionized). Add 8 drops from each of the indicator solutions. Titrate to pink. Note the amount in ml used ( $A_1$ ). Boil this solution for a few minutes. The solution will turn green. Cool rapidly to room temperature under running water. Continue the titration until the next pink colour change occurs. Note also this volume

( $A_2$ ). Boil the solution for a few minutes. Cool rapidly to room temperature under running water. Continue the titration until the next pink colour occurs. Note also this volume ( $A_3$ )

**Note!** Temperature changes will influence the volume and the concentration of the titrant solution. The working temperature of the titrant should approximate that of its temperature during standardization. If temperature corrections are necessary, sufficient accuracy may be obtained by use of a correction table. (AOAC 942.25)

### 7.2 Calculation

$$\text{Molarity (M)} = \frac{18,870 \times W_1}{(A_1 + A_2 + A_3)}$$

**Note!** Concentration must be accurate to four digits, i.e. 0.2000 M.

**Note!** The colour change of this official procedure (AOAC 936.15) may be difficult to see, therefore a pH meter or a mixed indicator (e.g. 0.1 g Methyl red and 0.1 g Bromocresol green in 100 ml methanol) will make it much easier to perform.

### 7.3 Receiver Solution

Boric acid 4 % with bromocresol green / methyl red indicator solution

In order to obtain accurate results the receiver solution is adjusted so that a small (0.05-0.20 ml) positive blank is obtained when running a blank sample. The 4 % boric acid receiver solution is prepared by dissolving 400 g of boric acid in about 5-6 l very hot deionized water. Mix and add more hot deionized water to a volume of about 9 l. Cool the solution to room temperature and add 100 ml of bromocresol green solution (100 mg in 100 ml methanol) and 70 ml of methyl red solution (100 mg in 100 ml of methanol). Dilute to 10 l with deionized water and mix carefully.

**Note!** The addition of alkali is to achieve a positive blank value. This should, however, be kept between 0.05 - 0.20 ml titrant, to obtain good repeatability when testing blanks.

Adjustment of the boric acid is made by the following procedure:

1. Transfer 25 ml boric acid solution to a receiver flask and add 100 ml of distilled water. If the solution in the flask is still red, titrate with 0.1 M sodium hydroxide solution until a neutral grey colour is obtained. Calculate the amount of sodium hydroxide solution necessary to adjust the boric acid solution in the 10 l flask with the formula: ml 1.0 M alkali = ml titrant x 40
2. Add the calculated amount of 1.0 M alkali solution to the boric acid solution. Mix.
3. To check proceed as follows using 25 ml of the boric acid solution. Run a blank. If the value of this blank is high (0.5 ml of 0.2 M HCl) the boric acid is incorrectly adjusted. This might create irregular blanks. For correction add HCl directly into the boric acid tank, mix it carefully and repeat until a reading of 0.05 - 0.20 ml HCl is obtained. If a positive blank is not achieved, add further small quantities of 1 M NaOH and repeat the check until a satisfactory value is achieved.

## Kjeltec™ 8100 Distillation Unit Tecator™ 2508/2520 Digestor

### 1 Scope

This PQ applies to the Digestion system 2508/2520 (including exhaust and scrubber unit) and Kjeltec 8100 Distillation Unit manufactured by FOSS Analytical. The user of the instrument performs the PQ.

### 2 Intended Use

The Digestion system (including exhaust and scrubber) and Kjeltec 8100 Distillation Unit are intended for laboratory use analyzing parameters as specified in FOSS Application Notes.

### 3 Purpose

The guidelines are intended to assist the user in successfully developing Performance Qualifications for the specific application(s) to which the instrument is applied.

The Performance Qualification (PQ) includes the process of demonstrating that the Digestion system 2508/2520 (including exhaust and scrubber unit) and the Kjeltec 8100 Distillation unit consistently perform according to a specification appropriate for its routine use. Main activities in the PQ phase are:

- Preventive maintenance
- On-going verification tests

This document suggests routines to fulfill the requirements for an acceptable PQ but the final procedure should be adapted to local routines for similar equipment.

### 4 Definition of Test Procedures

#### 4.1 Preventive Maintenance

Maintenance of the Kjeltec 8100 should be performed according to the instructions in manual, see User Manual Kjeltec 8100/8200 Distillation Unit, chapter 5. Maintenance. A yearly service is recommended (service agreement).

Maintenance of the Digestion block (including exhaust and scrubber) should be performed according to instruction in the user manual, see User Manual Tecator Digestor, chapter 5. Maintenance.

#### Dedicated Analytical Solutions

FOSS Analytical A/S  
69 Slangerupgade  
DK-3400 Hillerød  
Denmark

Tel +45 7010 3370  
Fax +45 7010 3371  
E-mail support@foss.dk  
Web www.foss.dk

FOSS Analytical AB  
Box 70  
SE-263 21 Höganäs  
Sweden

Tel +46 42 361500  
Fax +46 42 340349  
E-mail support@foss.dk  
Web www.foss.dk

### 4.2 Ongoing Qualification Tests

#### Block Temperature

The temperature for the digestion is limited by the boiling point for the sulphuric acid, this can be increased by adding a salt ( $K_2SO_4$ ) to the digestion mixture. It's important that the optimal ratio between acid and salt is kept; please follow recommendation in AN 300 or suggested procedures for a specific kind of sample material.

The block temperature itself can be controlled external by inserting a temperature probe in the intended hole in the aluminium block (front row of holes).

Use the reagents and method procedure specified in AN 300. Use only reagents of recognized analytical grade, unless otherwise specified and distilled or demineralised water or water of equivalent purity.

*Suggested standard material for internal quality control:*

Ammonium sulphate  $[(NH_4)_2SO_4]$ , min. 99.5 % (mass fraction), with certified purity.

**Note: The above chemical is usually readily available with a certificate specifying the purity.**

Alternatively ammonium iron(II) sulphate,  $(NH_4)_2Fe(SO_4)_2 \times 6 H_2O$ , with certified purity may be used.

Tryptophan ( $C_{11}H_{12}N_2O_2$ ), minimum assay 99 % (mass fraction). Nitrogen content 137.2 g/kg. Do not dry in an oven before use.

Acetanilide ( $C_8H_9NO$ ), minimum assay 99 % (mass fraction). Nitrogen content 103.6 g/kg. Do not dry in an oven before use.

Sucrose, ( $C_{12}H_{22}O_{11}$ ), with a nitrogen content of not more than 0.002 % (mass fraction). Do not dry in an oven before use.

#### Blank Tests

Carry out a blank test following the currently used procedure for digestion, distillation and titration taking 2 ml of water and about 0.7 g of sucrose instead of the test portion. Keep a record of blank values. If blank values change, identify the cause.

**Note: The amount of titrant used in the blank test should always be greater than 0.0 ml. Blanks within the same laboratory should be consistent across time.**

### 4.3 Recovery Tests

Regularly run recovery studies to check the accuracy of procedure and equipment:

- *Nitrogen loss.* - Use 0.12 g ammonium sulphate and 0.67 g sucrose per flask weighed to the nearest 0.1 mg. Add all other reagents as stated in the method currently used (Kjeltabs,  $H_2SO_4$ , etc). Digest and distil under same conditions as for sample. Recoveries shall be >99 %.
- *Digestion efficiency* - Use a test portion of minimum 0.15 g of tryptophan or acetanilide and 0.67 g sucrose per flask weighed to the nearest 0.1 mg. Determine the nitrogen content according to the current procedure in use. The recoveries of tryptophan shall be >98.5 %; the recoveries of acetanilide shall be >99.5 %.
- *Distillation and titration efficiency* - Distil 0.10 – 0.15 g  $\pm 0.0001$  g ammonium sulphate, omitting the digestion step. The recoveries should be >99.5 %.

**Note: Results less than 98.5 % or more than 101.0 % in either of the recovery tests indicate failures in the procedure and/or inaccurate concentration of the standard volumetric hydrochloric acid solution (should be adjusted to four decimals accuracy according to procedure in AN 300)**

External Quality Control Program

It is recommended to participate in an external quality control program, such a proficiency program or ring test, with equivalent sample material as analysed within the laboratory.

Calculation and Expression of Results

Equation for w\_n calculation: w\_n = (14.007(V\_s - V\_b)N x 100%) / m

- Where: w\_n is the nitrogen content of the sample, expressed as a percentage by mass.
- V\_s is the numerical value of the volume of the hydrochloric acid standard volumetric solution) used in the sample test, in milliliters, expressed to the nearest 0.05 ml.
- V\_b is the numerical value of the volume of the hydrochloric acid standard volumetric solution used in the blank test, in milliliters, expressed to the nearest 0.05 ml.
- N is the numerical value of the exact normality of the hydrochloric acid standard volumetric solution, expressed to four decimal places.
- m is the numerical value of the mass of the test portion , in milligrams, expressed to the nearest 1 mg for sample weights >1 g or to the nearest 0.1 mg for sample weights <1 g.

5 Maintenance

5.1 Maintenance Kjeltect™ 8100

See instructions in User Manual - Kjeltect 8100/8200, chapter 5 Maintenance.

5.2 Maintenance Tecator™ Digestor

See instructions in User Manual - Tecator Digestor, chapter 5 Maintenance.

6 The Maintenance Record Charts

This record charts are provided to assist you in keeping your system in good working order. Please make copies and use them regularly as they can often help us to help you in the unlikely event a system malfunction.

เอกสารไม่ควบคุม

6.1 FossCare™ Service Log

Table with 5 columns: Date, Work Performed and Parts Used, FOS Representative, Customer Representative, FossCare™ Service Report No. The table contains 15 empty rows for recording service events.

เอกสารไม่ควบคุม

Applicable for FOSS sales and service companies.



6.2.2 Weekly Maintenance

Date	/	/	/	/	/	/	/	/	/	/	/	/	/	/	/	/	/	/	/	/
Check of Safety Functions																				
Check of Reagent Tanks																				
Check of Digestion Tubes																				
Signature																				
Number of analyses																				

เอกสารไม่ควบคุม

6.2.3 Every 1-3 Months Maintenance

Date	/	/	/	/	/	/	/	/	/	/	/	/	/	/	/	/	/	/	/	/
Cleaning of Alkali Pump																				
Cleaning of Splash Head																				
Cleaning of Receiver Solution Dispensing System																				
Check of Tube Adapter																				
Signature																				
Number of analyses																				

เอกสารไม่ควบคุม



[illegible]

11(11)

## เอกสารไม่ควบคุม

DQE Services



NSC-TISI-TIS 17025  
CALIBRATION 0404

**เอกสารไม่ควบคุม**

FM-708-02 R01 1/11/2021

## REPORT OF CALIBRATION

Certificate No. : SP22-016

Page 2 of 5

Environment Condition : Ambient Temperature  $25 \pm 5$  °CRelative humidity  $55 \pm 20$  %RH

Calibration method : In-house method CP-01 Based on ASTM E275-08

## Certified Reference Materials :

Material	Serial No.	Certificate No.	Due date
Absorbance Standard set	25760	95935	22 October 2023
Absorbance Standard set	25757	95929	22 October 2023
Wavelength Standard set	25806	95916	22 October 2023
Wavelength Standard set	25758	95915	22 October 2023

Traceability : This certification is traceable to the International System of Unit maintained at National -

Institute of Standards and Technology (NIST) through Starna Scientific Limited

Spectral Band Width of UUC : 1.5 nm.

Scan Speed of UUC : 90 nm/min

Scan Interval of UUC : 0.15 nm.

Resolution of UUC : Photometric 0.0001 Abs.

Wavelength 0.1 nm.

เอกสารไม่ควบคุม

## REPORT OF CALIBRATION

Certificate No. : SP22-016

Page 3 of 5

Calibration Results : Without adjustment

## Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor <i>k</i>
420	0.0000	0.0000	0.0000	0.0028	2.00
	0.5787	0.5755	0.0032	0.0031	2.00
	1.0490	1.0436	0.0054	0.0029	2.00
	2.1900	2.1847	0.0053	0.0075	2.00
440	0.0000	0.0000	0.0000	0.0028	2.00
	0.5607	0.5588	0.0019	0.0034	2.00
	1.0247	1.0232	0.0015	0.0035	2.00
	2.1229	2.1211	0.0018	0.0082	2.00
465	0.0000	0.0000	0.0000	0.0028	2.00
	0.5236	0.5197	0.0039	0.0029	2.00
	0.9634	0.9625	0.0009	0.0028	2.00
	1.9763	1.9752	0.0011	0.0070	2.00
546.1	0.0000	-0.0001	0.0001	0.0028	2.00
	0.5191	0.5171	0.0020	0.0031	2.00
	1.0003	0.9984	0.0019	0.0033	2.00
	1.9987	1.9946	0.0041	0.0084	2.00
590	0.0000	0.0000	0.0000	0.0028	2.00
	0.5523	0.5509	0.0014	0.0030	2.00
	1.0809	1.0799	0.0010	0.0029	2.00
	2.0391	2.0329	0.0062	0.0080	2.00
635	0.0000	0.0000	0.0000	0.0028	2.00
	0.5601	0.5584	0.0017	0.0031	2.00
	1.0512	1.0498	0.0014	0.0029	2.00
	1.9294	1.9265	0.0029	0.0082	2.00

เอกสารไม่ควบคุม

**REPORT OF CALIBRATION**

Certificate No. : SP22-016

Page 4 of 5

**Photometric Accuracy :**

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor <i>k</i>
235	0.0000	0.0001	-0.0001	0.0050	2.00
	0.7478	0.7421	0.0057	0.0056	2.00
257	0.0000	0.0000	0.0000	0.0050	2.00
	0.8686	0.8619	0.0067	0.0059	2.00
313	0.0000	0.0000	0.0000	0.0050	2.00
	0.2912	0.2896	0.0016	0.0051	2.00
350	0.0000	0.0000	0.0000	0.0050	2.00
	0.6448	0.6403	0.0045	0.0055	2.00

**เอกสารไม่ควบคุม**

FM-708-02 R01 1/11/2021

**REPORT OF CALIBRATION**

Certificate No. : SP22-016

Page 5 of 5

**Wavelength Accuracy :**

CRMs Values (nm.)	UUC Reading (nm.)	Correction (nm.)	Uncertainty (nm.)	Coverage factor <i>k</i>
241.72	242.0	-0.28	0.18	2.00
279.45	279.5	-0.05	0.18	2.00
287.81	287.5	0.31	0.18	2.00
334.06	333.5	0.56	0.18	2.00
360.93	360.5	0.43	0.18	2.00
418.59	418.0	0.59	0.18	2.00
445.94	445.4	0.54	0.18	2.00
453.66	453.2	0.46	0.18	2.00
460.02	459.7	0.32	0.18	2.00
536.59	536.2	0.39	0.18	2.00
637.98	638.3	-0.32	0.18	2.00
431.38	431.0	0.38	0.18	2.00
472.50	472.5	0.00	0.18	2.00
513.47	513.5	-0.03	0.18	2.00
528.88	528.5	0.38	0.18	2.00
573.17	573.0	0.17	0.18	2.00
585.35	585.0	0.35	0.20	2.00
684.40	684.7	-0.30	0.18	2.00
740.72	740.8	-0.08	0.20	2.00
748.55	748.5	0.05	0.18	2.00
807.03	807.3	-0.27	0.18	2.00
879.28	879.0	0.28	0.18	2.00

Remark : - UUC = Unit Under Calibration

- N/A = Not Available

- The result expanded uncertainty of measurement U is stated as the standard uncertainty of measurement multiplied by the coverage factor *k*,

which for a normal distribution corresponds to a coverage probability of approximately 95%

- \* Indicates non TISI accredited

- End of Certificate -

**เอกสารไม่ควบคุม**

FM-708-02 R01 1/11/2021

## CERTIFICATE OF CALIBRATION

Certificate No. : SP22-007

Page 1 of 5

Customer : United Analyst and Engineering Consultant Co.,Ltd. (Head Office)

Address : 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong,

Bangkok 10260

Location of calibration : Laboratory 315

Equipment : UV-Vis Spectrophotometer

Manufacturer : Hitachi

Model : U-1900

Serial No. : 2021-064

ID No. : UAE.WAS.006/2552


Received Date : 20 January 2022

Calibration Date : 20 January 2022

Issue Date : 24 January 2022

Condition Instrument : Good


Calibrated by :

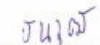

(Mr. Tanawut Rittidach)

Technical Manager

Approved by :

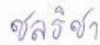

(Ms. Chonthicha Sangngern)

Quality Manager

The calibration result is applied only to the above calibrated item and was found accurate as shown on date and place of calibration only.

The measurement capability of the laboratory and its traceability to recognized national standards and to the unit of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the DQE Services Co., Ltd.

## REPORT OF CALIBRATION

Certificate No. : SP22-007

Page 2 of 5

Environment Condition : Ambient Temperature  $25 \pm 5$  °CRelative humidity  $55 \pm 20$  %RH

Calibration method : In-house method CP-01 Based on ASTM E275-08

Certified Reference Materials :

Material	Serial No.	Certificate No.	Due date
Absorbance Standard set	25760	95935	22 October 2023
Absorbance Standard set	25757	95929	22 October 2023
Wavelength Standard set	25806	95916	22 October 2023
Wavelength Standard set	25758	95915	22 October 2023

Traceability This certification is traceable to the International System of Unit maintained at National -

Institute of Standards and Technology (NIST) through Sarna Scientific Limited

Spectral Band Width of UUC : 4.0 nm.

Scan Speed of UUC : 200 nm/min

Scan Interval of UUC : 0.1 nm.

Resolution of UUC : Photometric 0.001 Abs.

Wavelength 0.1 nm.

**REPORT OF CALIBRATION**

Certificate No. : SP22-007

Page 3 of 5

Calibration Results : Without adjustment

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor <i>k</i>
420	0.0000	0.000	0.0000	0.0028	2.00
	0.5787	0.577	0.0017	0.0031	2.00
	1.0490	1.050	-0.0010	0.0029	2.00
	2.1900	2.183	0.0070	0.0080	2.00
440	0.0000	0.000	0.0000	0.0028	2.00
	0.5607	0.560	0.0007	0.0034	2.00
	1.0247	1.023	0.0017	0.0035	2.00
	2.1229	2.118	0.0049	0.0079	2.00
465	0.0000	0.000	0.0000	0.0028	2.00
	0.5236	0.521	0.0026	0.0030	2.00
	0.9634	0.963	0.0004	0.0029	2.00
	1.9763	1.974	0.0023	0.0070	2.00
546.1	0.0000	0.000	0.0000	0.0028	2.00
	0.5191	0.518	0.0011	0.0031	2.00
	1.0003	1.000	0.0003	0.0033	2.00
	1.9987	1.996	0.0027	0.0084	2.00
590	0.0000	0.000	0.0000	0.0028	2.00
	0.5523	0.552	0.0003	0.0030	2.00
	1.0809	1.082	-0.0011	0.0030	2.00
	2.0391	2.033	0.0061	0.0079	2.00
635	0.0000	0.000	0.0000	0.0028	2.00
	0.5601	0.562	-0.0019	0.0031	2.00
	1.0512	1.052	-0.0008	0.0030	2.00
	1.9294	1.925	0.0044	0.0079	2.00

เอกสารไม่ควบคุม

**REPORT OF CALIBRATION**

Certificate No. : SP22-007

Page 4 of 5

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor <i>k</i>
235	0.0000	0.000	0.0000	0.0050	2.00
	0.7478	0.746	0.0018	0.0057	2.00
257	0.0000	0.000	0.0000	0.0050	2.00
	0.8686	0.861	0.0076	0.0059	2.00
313	0.0000	0.000	0.0000	0.0050	2.00
	0.2912	0.291	0.0002	0.0051	2.00
350	0.0000	0.000	0.0000	0.0050	2.00
	0.6448	0.638	0.0068	0.0055	2.00

เอกสารไม่ควบคุม

## REPORT OF CALIBRATION

Certificate No. : SP22-007

Page 5 of 5

## Wavelength Accuracy :

CRMs Values (nm.)	UUC Reading (nm.)	Correction (nm.)	Uncertainty (nm.)	Coverage factor k
241.54	240.8	0.74	0.18	2.00
279.40	278.5	0.90	0.18	2.00
288.70	288.0	0.70	0.18	2.00
334.22	333.5	0.72	0.18	2.00
361.26	360.5	0.76	0.18	2.00
418.48	418.0	0.48	0.18	2.00
446.70	446.0	0.70	0.18	2.00
453.20	453.0	0.20	0.18	2.00
460.06	459.5	0.56	0.18	2.00
536.90	536.0	0.90	0.18	2.00
637.94	637.2	0.74	0.18	2.00
440.74	440.0	0.74	0.18	2.00
472.22	471.6	0.62	0.18	2.00
513.70	513.0	0.70	0.18	2.00
528.72	528.0	0.72	0.18	2.00
574.60	573.8	0.80	0.18	2.00
585.48	584.6	0.88	0.20	2.00
684.63	684.0	0.63	0.18	2.00
740.27	739.8	0.47	0.20	2.00
748.28	747.8	0.48	0.18	2.00
807.16	806.4	0.76	0.18	2.00
879.70	878.8	0.90	0.18	2.00

Remark : - UUC = Unit Under Calibration

- N/A = Not Available

- The result expanded uncertainty of measurement U is stated as the standard uncertainty of measurement multiplied by the coverage factor k, which for a normal distribution corresponds to a coverage probability of approximately 95%

- \* Indicates non TISI accredited

- End of Certificate -

FM-708-02 R01 1/11/2021

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
 CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES  
 534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250  
 TEL. 0-2717-3000-27 FAX. 0-2719-9484



Cert. No.: 22TM563

Page.: 1 of 3

## Certificate of Calibration

Equipment : Incubator

Manufacturer : Memmert

Model : IPP 260

Serial No. : V615.0187

ID No. : UAE.MIC.003/2559

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.  
 3 Soi Udomsuk 41, Sukhumvit Road,  
 Bangchak, Phrakhanong,  
 Bangkok 10260

Location : Microbiology Laboratory

Received Order : 7 April 2022

Calibration Date : 7 April 2022

Ambient Temperature : ( 26 ± 10 ) °C

Relative Humidity : ( 50 ± 30 ) %

Calibrated by : Prawit Sodavitchit

Approved by :

( ) Pornthippa Tameyakul  
 ( ) Malee Butkruea  
 ( ) Suwit Imjai

Issue Date : 18 April 2022

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
 Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม

A 0040248



Equipment : Incubator  
 Condition As-Received : Used Item  
 Reference : 2204-0016OC-1

Cert. No.: 22TM563  
 Page.: 2 of 3

**Procedure Used :-**

Calibration were conducted using calibration procedure CP-OT02 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector ( RTD ).

The temperature scale used was based on ITS-90.

**Condition of this result of calibration**

1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1 ) Data Acquisition	34970A	MY44067817	21LM10	20 Jul 2022


2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

**Result of Calibration :-** ( \* ) Without Adjustment

**Function of UUC\* :** Temperature Source

**Fresh air setting :** Close

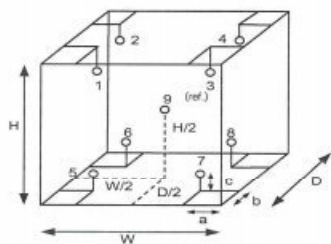

**Probe Installation Details :**

a = 5.0 cm  
 b = 5.0 cm  
 c = 5.0 cm

**Dimension of Chamber :**

D = 0.50 m  
 W = 0.64 m  
 H = 0.80 m  
 Capacity = 0.26 m<sup>3</sup>

Environment during calibration		
	Beginning	Finished
Temp. ( °C )	26	26
REL.Humid. ( % )	60	62
AC Supply ( Volt )	220	220

Position :	Ref. Std. ID No.:
1	15RTD2/11
2	15RTD2/12
3	15RTD2/13
4	15RTD2/14
5	15RTD2/15
6	15RTD2/16
7	15RTD2/17
8	15RTD2/18
9 (ref.)	15RTD2/19



Equipment : Incubator  
 Condition As-Received : Used Item  
 Reference : 2204-0016OC-1

Cert. No.: 22TM563  
 Page.: 3 of 3

**Result of Calibration :-** ( \* ) Without Adjustment

**Function of UUC\* :** Temperature Source

**Fresh air setting :** Close

Calibration Point ( °C )	UUC* Setting ( °C )	UUC* Reading ( °C )	Temperature stability ( ± °C )	Temperature uniformity ( °C )	Overall Variation ( °C )	Uncertainty ( ± °C )	Coverage Factor k
35.0	35.0	35.0	0.12	0.53	0.79	0.30	2

Calibration Point ( °C )	Measured Temperature ( °C )								
	Position								
	1	2	3	4	5	6	7	8	9 (ref.)
35.0	35.170	35.167	34.938	34.844	34.816	34.854	34.584	34.730	34.780

**Average\* :** The average of 30 values in each position.

**Temperature stability :** One-half of the greatest maximum difference of measured temperature at any one sensor.

**Temperature uniformity :** The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

**Overall Variation :** The Difference of the maximum and minimum measured temperatures throughout observation.

**UUC\* :** Unit Under Calibration

**Note :** The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่ควบคุม

a 1104310

เอกสารไม่ควบคุม

a 1104309



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES  
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250  
TEL. 0-2717-3000-27 FAX. 0-2719-9484



Cert. No.: 22TM503  
Page.: 1 of 3

## Certificate of Calibration

**Equipment :** Incubator  
**Manufacturer :** Memmert  
**Model :** IPP 260  
**Serial No. :** V618.0033  
**ID No. :** UAE.MIC.021/2561  
**Submitted by :** United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260  
**Location :** Microbiology Laboratory (302)  
**Received Order :** 3 May 2022  
**Calibration Date :** 3 May 2022  
**Ambient Temperature :**  $(26 \pm 10) ^\circ\text{C}$   
**Relative Humidity :**  $(50 \pm 30) \%$

**Calibrated by :** Krisda Malee

**Approved by :**

Approved Signatory

( ) Pornthippa Tameyakul  
(x) Malee Butkruea  
( ) Suwit Imjai

**Issue Date :** 10 May 2022

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3: Equipment Calibration and Testing

เอกสารไม่ควบคุม



**Equipment :** Incubator  
**Condition As-Received :** Used Item  
**Reference :** 2205-0003OC-4  
**Procedure Used :-**

**Cert. No.:** 22TM503  
**Page.:** 2 of 3

Calibration were conducted using calibration procedure CP-OT02 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector ( RTD ).

The temperature scale used was based on ITS-90.

### Condition of this result of calibration

#### 1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1 ) Data Acquisition	34970A	MY49023932	21LM8	06 Jul 2022


2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

**Result of Calibration :-** ( \* ) Without Adjustment

**Function of UUC\* :** Temperature Source

**Fresh air setting :** Not Available

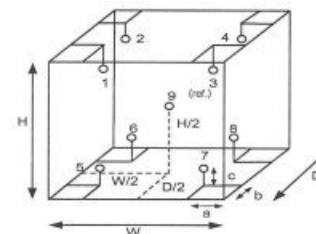

#### Probe Installation Details :

a = 5.0 cm  
b = 5.0 cm  
c = 5.0 cm

#### Dimension of Chamber :

D = 0.50 m  
W = 0.64 m  
H = 0.80 m  
Capacity = 0.26 m<sup>3</sup>

Environment during calibration		
	Beginning	Finished
Temp. ( °C )	23	22
REL.Humid. ( % )	53	54
AC Supply ( Volt )	221	220

Position :	Ref. Std. ID No.:
1	20-16RTD-01
2	20-16RTD-02
3	20-16RTD-03
4	20-16RTD-04
5	20-16RTD-05
6	20-16RTD-06
7	20-16RTD-07
8	20-16RTD-08
9 (ref.)	20-16RTD-09

เอกสารไม่ควบคุม



Equipment : Incubator  
Condition As-Received : Used Item  
Reference : 2205-0003OC-4  
Result of Calibration :- ( \* ) Without Adjustment  
Function of UUC\* : Temperature Source  
Fresh air setting : Not Available

Cert. No.: 22TM503  
Page.: 3 of 3

Calibration Point ( °C )	UUC* Setting ( °C )	UUC* Reading ( °C )	Temperature stability ( ± °C )	Temperature uniformity ( °C )	Overall Variation ( °C )	Uncertainty ( ± °C )	Coverage Factor <i>k</i>
22.0	22.0	22.0	0.051	0.095	0.19	0.30	2
44.0	44.0	44.0	0.10	0.83	1.2	0.32	2

Calibration Point ( °C )	Measured Temperature ( °C )								
	Position								
	1	2	3	4	5	6	7	8	9 (ref.)
22.0	22.238	22.200	22.139	22.207	22.184	22.178	22.136	22.169	22.161
44.0	44.604	44.670	44.240	44.315	43.974	44.446	43.584	44.209	44.306

Average\* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC\* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES  
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250  
TEL. 0-2717-3000-27 FAX. 0-2719-9484



Cert. No.: 22TM565  
Page.: 1 of 3

## Certificate of Calibration

Equipment : Water Bath  
Manufacturer : Memmert  
Model : WNE 14  
Serial No. : L414.1407  
ID No. : UAE.MIC.006/2558  
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260  
Location : Microbiology Laboratory  
Received Order : 7 April 2022  
Calibration Date : 7 April 2022  
Ambient Temperature : ( 26 ± 10 ) °C  
Relative Humidity : ( 50 ± 30 ) %

Calibrated by :

Approved by :

( ) Pornthippa Tameyakul  
( ✓ ) Malee Butkruea  
( ) Suwit Imjai

Issue Date : 18 April 2022

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing

เอกสารไม่ควบคุม



Equipment : Water Bath  
Condition As-Received : Used Item  
Reference : 2204-0016OC-4  
Procedure Used :-

Cert. No.: 22TM565  
Page.: 2 of 3

Calibration were conducted using in-house calibration procedure CP-OT04 according to direct measurement method with Data Acquisition which connected with Industrial Platinum Resistance Thermometer ( IPRT ).

The temperature scale used was based on ITS-90.

**Condition of this result of calibration**

1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1 ) Data Acquisition	34970A	MY44067817	21LM10	20 Jul 2022


2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

**Result of Calibration :-** ( \* ) Without Adjustment

**Function of UUC\* :** Temperature Source

	Environmental		AC Voltage Supply
	( °C )	( %R.H. )	( Volt )
Beginning of Calibration	26	62	220
Finished of Calibration	26	65	220

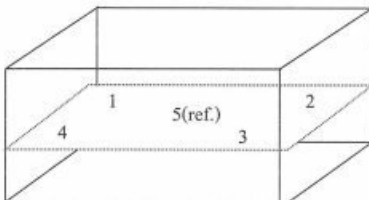

Front

Position :	Ref. Std. ID No.:
1	70RC143
2	70RC144
3	70RC145
4	70RC146
5(ref.)	70RC147

เอกสารไม่ควบคุม



Equipment : Water Bath  
Condition As-Received : Used Item  
Reference : 2204-0016OC-4  
Result of Calibration :- ( \* ) Without Adjustment  
Function of UUC\* : Temperature Source

Cert. No.: 22TM565  
Page.: 3 of 3

Calibration point ( °C )	UUC* Setting ( °C )	UUC* Reading ( °C )	Average* Standard Reading ( °C )				
			Position				
			1	2	3	4	5 (ref.)
44.5	44.5	44.5	44.424	44.409	44.478	44.470	44.581

Calibration point ( °C )	Uniformity ( °C )	Stability ( ± °C )	Uncertainty ( ± °C )	Coverage Factor k
44.5	0.22	0.039	0.15	2

**Average\* :** The average of 30 values in each position.

**Uniformity :** The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

**Stability :** One-half of the greatest maximum difference of measured temperature at any one probe.

**UUC\* :** Unit Under Calibration

**Note :** The reported uncertainty of measurement was included stability and excluded uniformity.

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES  
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250  
TEL. 0-2717-3000-27 FAX. 0-2719-9484



Cert. No.: 22TM564  
Page.: 1 of 3

## Certificate of Calibration

**Equipment :** Water Bath  
**Manufacturer :** Memmert  
**Model :** WNE 14  
**Serial No. :** L414.1410  
**ID No. :** UAE.MIC.007/2558  
**Submitted by :** United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260  
**Location :** Microbiology Laboratory  
**Received Order :** 7 April 2022  
**Calibration Date :** 7 April 2022  
**Ambient Temperature :** ( 26 ± 10 ) °C  
**Relative Humidity :** ( 50 ± 30 ) %

**Calibrated by :** Prawit Sodavitchit

**Approved by :**

( ✓ ) Pornthippa Tameyakul  
( ✓ ) Malee Butkruea  
( ) Suwit Imjai

**Issue Date :** 18 April 2022

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing

เอกสารไม่ควบคุม



**Equipment :** Water Bath  
**Condition As-Received :** Used Item  
**Reference :** 2204-0016OC-5  
**Procedure Used :-**

**Cert. No.:** 22TM564  
**Page.:** 2 of 3

Calibration were conducted using in-house calibration procedure CP-OT04 according to direct measurement method with Data Acquisition which connected with Industrial Platinum Resistance Thermometer ( IPRT ).

The temperature scale used was based on ITS-90.

### Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1 ) Data Acquisition	34970A	MY44067817	21LM10	20 Jul 2022


2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

**Result of Calibration :-** ( \* ) Without Adjustment

**Function of UUC\* :** Temperature Source

	Environmental		AC Voltage Supply
	( °C )	( %R.H. )	( Volt )
Beginning of Calibration	26	60	220
Finished of Calibration	26	62	220

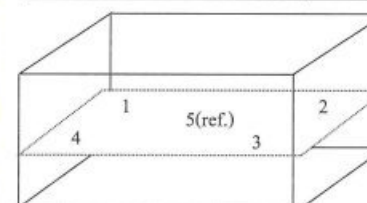

Front

Position :	Ref. Std. ID No.:
1	70RC143
2	70RC144
3	70RC145
4	70RC146
5(ref.)	70RC147

เอกสารไม่ควบคุม



Equipment : Water Bath  
Condition As-Received : Used Item  
Reference : 2204-0016OC-5  
Result of Calibration :- ( \* ) Without Adjustment  
Function of UUC\* : Temperature Source

Cert. No.: 22TM564  
Page.: 3 of 3

Calibration point ( °C )	UUC* Setting ( °C )	UUC* Reading ( °C )	Average* Standard Reading ( °C )				
			Position				
			1	2	3	4	5 (ref.)
44.5	44.5	44.5	44.498	44.530	44.542	44.635	44.591

Calibration point ( °C )	Uniformity ( °C )	Stability ( ± °C )	Uncertainty ( ± °C )	Coverage Factor <i>k</i>
44.5	0.16	0.068	0.15	2

**Average\*** : The average of 30 values in each position.

**Uniformity** : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

**Stability** : One-half of the greatest maximum difference of measured temperature at any one probe.

**UUC\*** : Unit Under Calibration

**Note** : The reported uncertainty of measurement was included stability and excluded uniformity.

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

เอกสาร

Calibration Certificate ID  
TH2058-096-040722-ACC-TH

METTLER TOLEDO

Mettler-Toledo (Thailand) Ltd.  
846/4 - 846/5 Lasalle Rd., Bangna Tai Sub-District  
Bangna District, Bangkok 10260  
+66 2723 0382  
MT-TH.ServiceSupport@mt.com



## Accuracy Calibration Certificate

### Customer

Company: United Analyst and Engineering Consultant Co., Ltd.  
Address: 3 Soi Udom Suk 41, Sukhumvit Rd., Bang Chak  
City: Phra Khanong Contact: Suwit Chotnok  
Zip / Postal: 10260  
State / Province: Bangkok  
Order Number: 

### Weighing Device

Manufacturer: Mettler Toledo Instrument Type: Weighing Instrument  
Model: MS603S/01 Asset Number: UAE.MIC.008/2553  
Serial No.: B007010311 Terminal Model: N/A  
Building: N/A Terminal Serial No.: N/A  
Floor: 2 Terminal Asset No.: N/A  
Room: Balance Room (206)

Range	Max. Capacity	Readability (d)
1	620 g	0.001 g

### Procedure

Calibration Guideline: EURAMET cg-18 v. 4.0 (11/2015)  
METTLER TOLEDO Work Instruction: CPW002/20

This calibration certificate contains measurements for As Found calibration. No As Left calibration was performed because the device was not modified after As Found calibration. Therefore, results for As Left correspond to As Found.

The sensitivity/span of the weighing instrument was adjusted before calibration with a built-in weight.

In accordance with EURAMET cg-18 (11/2015), the test loads were selected to reflect the specific use of the weighing device or to accommodate specific calibration conditions.

	Temperature		Humidity	
As Found	Start: 22.8 °C	End: 23.0 °C	Start: 49.9 %	End: 58.3 %

As Found Calibration Date: 07-Apr-2022 Calibrator:  
As Left Calibration Date: N/A  
Issue Date: 08-Apr-2022  
Approved Signatory:

เอกสารไม่ควบคุม


## Measurement Results

### Repeatability

Test Load: 200 g

	As Found	As Left
1	200.001 g	N/A
2	200.001 g	N/A
3	200.001 g	N/A
4	200.001 g	N/A
5	200.001 g	N/A
6	200.000 g	N/A
7	200.001 g	N/A
8	200.001 g	N/A
9	200.000 g	N/A
10	200.001 g	N/A

Standard Deviation	0.0004 g	N/A
--------------------	----------	-----

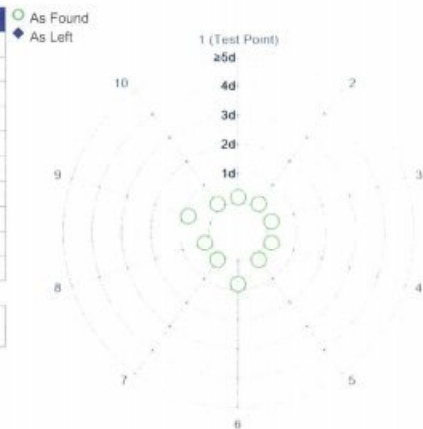

The "d" in the graph represents the readability of the range/interval in which the test was performed.


The results of this graph are based upon the absolute values of the differences from the mean value.

### Eccentricity

Test Load: 200 g

Position	As Found	As Left
1	200.001 g	N/A
2	200.001 g	N/A
3	200.002 g	N/A
4	200.002 g	N/A
5	200.000 g	N/A

Maximum Deviation	0.001 g	N/A
-------------------	---------	-----

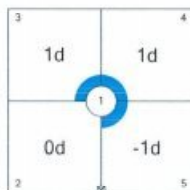


As Found

The "d" in the graph represents the readability of the range/interval in which the test was performed.

## Error of Indication

As Found

	Reference Value	Indication	Error of Indication	Expanded Uncertainty	k
1	0.000 g	0.000 g	0.000 g	1.0 mg	2
2	0.500 g	0.500 g	0.000 g	1.2 mg	2
3	1.000 g	1.000 g	0.000 g	1.2 mg	2
4	50.000 g	50.000 g	0.000 g	1.2 mg	2
5	100.000 g	100.000 g	0.000 g	1.3 mg	2
6	150.000 g	150.000 g	0.000 g	1.5 mg	2
7	200.000 g	200.001 g	0.001 g	1.6 mg	2
8	300.001 g	300.001 g	0.000 g	2.0 mg	2
9	400.001 g	400.001 g	0.000 g	2.5 mg	2
10	500.001 g	500.002 g	0.001 g	2.9 mg	2
11	600.001 g	600.001 g	0.000 g	3.4 mg	2

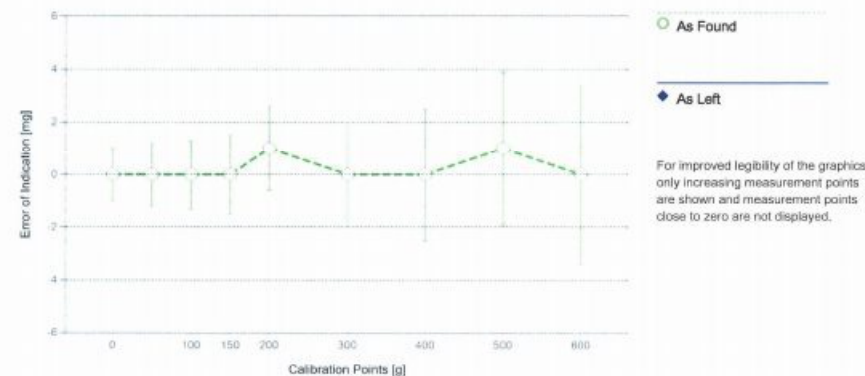

The uncertainty stated is the expanded uncertainty at calibration obtained by multiplying the standard combined uncertainty by the coverage factor  $k = 2$  which can be larger than 2 according to EURAMET cg-18. The value of the measurand lies within the assigned range of values with a probability of approximately 95%.

The user is responsible for maintaining environmental conditions and the settings of the weighing instrument when it was calibrated.

### Test Equipment

All weights used for metrological testing are traceable to national or international standards. The weights were calibrated and certified by an accredited calibration laboratory.

#### Weight Set 1: OIML F1

Weight Set No.: WS55 Date of Issue: 09-Jul-2021  
Certificate Number: CCM-0137-21-C Calibration Due Date: 07-Jul-2022

#### Weight Set 2: OIML E2

Weight Set No.: WS80 Date of Issue: 23-Feb-2022  
Certificate Number: C208581831 Calibration Due Date: 14-Aug-2023

#### Thermo Hygrometer

Equipment No.: IN161 Date of Issue: 14-Jun-2021  
Certificate Number: 21H1220 Calibration Due Date: 01-Jun-2022

### Remarks

FACT adjustment functionality activated  
Equipment condition: Good  
Next calibration according to customer's procedure  
Calibration data not decide by calibration laboratory

#### End of Accredited Section

The information below and any attachments to this calibration certificate are not part of the accredited calibration.

### Measurement Uncertainty of the Weighing Instrument in Use

Stated is the expanded uncertainty with  $k=2$  in use. The formula shall be used for the estimation of the uncertainty under consideration of the errors of indication. The value  $R$  represents the net load indication in the unit of measure of the device.

Temperature coefficient for the evaluation of the measurement uncertainty in use:  $3.0 \cdot 10^{-6} / K$

Temperature range on site for the evaluation of the measurement uncertainty in use: 3 K



#### Linearization of Uncertainty Equation

Range	As Found		As Left
	d	Max	
1	0.001 g	620 g	$U_1 = 1.2 \text{ mg} + 0.0186 \text{ mg/g} \cdot R$
			N/A

To optimize the stability of the linearization, besides of the zero load only increasing measurement points with a test load of 5% of the measurement range or larger are taken for the calculation of the linear equation.

#### Absolute and Relative Measurement Uncertainty in Use for Various Net Indications (Examples)

Net Indication	As Found		As Left	
0.062 g	1.2 mg	1.9%	N/A	N/A
0.620 g	1.2 mg	0.20%	N/A	N/A
6.200 g	1.3 mg	0.021%	N/A	N/A
62.000 g	2.4 mg	0.0038%	N/A	N/A
620.000 g	13 mg	0.0021%	N/A	N/A

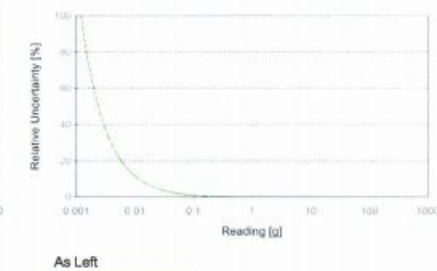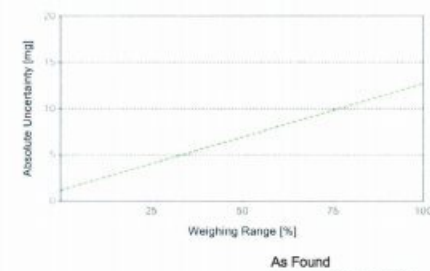



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES  
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250  
TEL. 0-2717-3000-27 FAX. 0-2719-9484



## Certificate of Calibration

Cert. No.: 22TM89  
Page.: 1 of 3

Equipment : Autoclave  
Manufacturer : ALP  
Model : CL-40L  
Serial No. : 802664  
ID No. : UAE.MIC.014/2550  
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260  
Location : Air Analysis Unit  
Received Order : 17 February 2022  
Calibration Date : 17 February 2022  
Ambient Temperature : ( 26 ± 10 ) °C  
Relative Humidity : ( 50 ± 30 ) %  
Calibrated by : Kunchit Promprat  
Approved by :   
( ) Pornthippa Tameyakul  
( ) Malee Butkruea  
( ) Suwit Imjai  
Issue Date : 22 February 2022

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3: Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : Autoclave  
Condition As-Received : Used Item  
Reference : 2202-0444OC-1  
Procedure Used :-

Cert. No.: 22TM89  
Page.: 2 of 3

Calibration were conducted using in-house calibration procedure CP-OT03 according to direct measurement method with Data Acquisition which connected with Thermocouple Type T

The temperature scale used was based on ITS-90.

### Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1 ) Data Acquisition	34970A	MY44035217	21LM30	23 Dec 2022

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

4. This result of calibration covers laboratory autoclaves for the sterilization of goods and material which could be infected with organisms categorized as Hazard Group 1, 2 and 3\*\*


(\*\* = Categorization of pathogens according to hazard and categories of containment, second edition, 1990 )

It does not cover autoclaves for use with material infect with organisms in Hazard Group 4, for which complete containment and sterilization of infected condensate is considered to be essential.

This result of calibration does not apply to sterilizers or disinfectors used for medical, dental, pharmaceutical or veterinary purposes which are directly concerned with patient care, or those used for fabrics subjected to sterilization which are required to be dry at the end of cycle.

Result of Calibration :- ( \* ) Without Adjustment

Function of UUC\* : Temperature Source

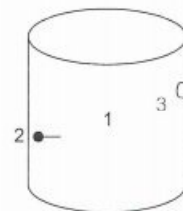

	Environmental		
	( °C )	( %R.H. )	( Volt )
Beginning of Calibration	27	68	226
Finished of Calibration	27	65	226

Position	Description	Ref. Std. ID No.:
1 =	Center of chamber	22-10TC-01
2 =	Temperature sensor	22-10TC-02
3 =	Exhaust port	22-10TC-03

เอกสารไม่ควบคุม



Equipment : Autoclave

Condition As-Received : Used Item

Reference : 2202-0444OC-1

Result of Calibration :- ( \* ) Without Adjustment

Cert. No.: 22TM89

Page.: 3 of 3

Operating parameter Set : Temperature = 122 °C

Sterilization period = 30 minute

UUC* Setting ( °C )	UUC* Reading ( °C )	Position	Average* Standard Reading ( °C )	Stability ( ± °C )	Pressure Reading ( MPa )	Uncertainty ( ± °C )	Coverage Factor <i>k</i>
122	122	1	122.373	0.32	0.12	1.2	2
		2	122.421				
		3	122.292				

Average\* : The average of 30 values in each position.

Stability : One-half of the greatest maximum difference of measured temperature at any one probe.

UUC\* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

เอกสาร



ภาคผนวก 13      เอกสารการขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

---





ที่ อก ๐๓๑๐(๑)/ ๑๘๗ ๙

กรมโรงงานอุตสาหกรรม  
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท  
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๐ ๙ กุมภาพันธ์ ๒๕๖๕

เรื่อง ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

เรียน กรรมการผู้จัดการ บริษัท ยูโนเต็ด แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด

อ้างถึง คำขอขึ้นทะเบียน/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน  
ลงวันที่ ๒๗ ธันวาคม ๒๕๖๔

สิ่งที่ส่งมาด้วย ๑. รายชื่อผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๔๐ ราย  
๒. รายชื่อเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๑๐๖ ราย  
๓. ขอบข่ายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม

ตามหนังสือที่อ้างถึง บริษัท ยูโนเต็ด แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด ขอต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ๖-๑๔๕ สถานที่ตั้งเลขที่ ๓ ซอยอุดมสุข ๔๔ ถนนสุขุมวิท แขวงบางจาก เขตพระโขนง กรุงเทพมหานคร ต่อกรมโรงงานอุตสาหกรรม นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้บริษัท ยูโนเต็ด แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน โดยมีองค์ประกอบดังนี้  
ก. ผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๔๐ ราย ตามสิ่งที่ส่งมาด้วย ๑  
ข. เจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๑๐๖ ราย ตามสิ่งที่ส่งมาด้วย ๒  
ค. ขอบข่ายสารมลพิษที่ได้รับขึ้นทะเบียนให้วิเคราะห์ในน้ำเสีย น้ำใต้ดิน อากาศเสีย สิ่งปฏิกูล หรือวัสดุที่ไม่ใช่แล้ว และดิน ตามสิ่งที่ส่งมาด้วย ๓

หนังสือฉบับนี้จะหมดอายุในวันที่ ๒ กุมภาพันธ์ ๒๕๖๘ หากประสงค์จะต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ให้ยื่นคำขอต่ออายุพร้อมเอกสารประกอบคำขอต่อกกรมโรงงานอุตสาหกรรมภายใน ๓๐ วัน ก่อนวันสิ้นสุดของหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ทั้งนี้ สามารถยื่นคำขอผ่านระบบอิเล็กทรอนิกส์ได้ที่หน้าเว็บไซต์กรมโรงงานอุตสาหกรรม ตาม QR Code ท้ายหนังสือฉบับนี้

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ



ผู้อำนวยการกองวิจัยและเตือนภัยมลพิษโรงงาน  
ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม

ยื่นคำขอผ่านระบบอิเล็กทรอนิกส์

กองวิจัยและเตือนภัยมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕

โทรสาร ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๕๙

ไปรษณีย์อิเล็กทรอนิกส์ saraban@diw.gmail.go.th

สิ่งที่ส่งมาด้วย ๑

เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท ยูโนเต็ด แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด เลขทะเบียน ๖-๑๔๕

ที่ อก ๐๓๑๐(๑)/ ๑๘๗ ๙

ลงวันที่ ๐ ๙ กุมภาพันธ์ ๒๕๖๕

ก. ผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๔๐ ราย

๑) นางสาวกฤษวรรณ กัทธิธิกุล

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๐๑

๒) นายณรงค์ ติมพาลี

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๐๒

๓) นางสาวนันติดา บุญไชย

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๐๓

๔) นางปิยะพัชร สุทมนัสวงษ์

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๐๔

๕) นางมานิดา แยมโย

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๐๕

๖) นางสาวเบญจวรรณ วิริยทัย

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๐๖

๗) นายพนรัตน์ วงศ์อนุรักษชัย

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๐๗

๘) นางสาวฉวีวรรณ บุญลา

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๐๘

๙) นายสุวิทย์ จอดนอก

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๐๙

๑๐) นางสาวโชติกา สมบูรณ์

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๑๐

๑๑) นางสาวบุษกร เลิศภาณุมาศ

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๑๑

๑๒) นางสาววิไลลักษณ์ ศรีสุข

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๑๒

๑๓) นางสาวปิรญา จรัสโชติพันธ์

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๑๓

๑๔) นายศิลา บรรจงใจรักษ์

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๑๔

๑๕) นายปฏิกรณ์ คณะนา

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๑๕

๑๖) นายธีรวัฒน์ ชมมิ่ง

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๑๖

๑๗) นางสาวศิริพร ศรีประดิษฐ์

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๑๗

๑๘) นางสาวสาวตรี วีรัง

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๑๘

๑๙) นางสาวนพวรรณ อูราภักษ์

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๑๙

๒๐) นายภูษงค์ พานิชย์เลิศอำไพ

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๒๐

๒๑) นายณัฐวัฒน์ แดงสวัสดิ์

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๒๑

๒๒) นายเอกรัตน์ ปณะคามินทร์

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๒๒

๒๓) นางสาวนิศารัตน์ ศรีสกุลสิทธิโชค

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๒๓

๒๔) นางสาวเจตจรินทร์ ทำสะอาด

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๒๔

๒๕) นางสาวสุวรรณ คงทอง

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๒๕

๒๖) นางสาววรกร พัดสองชั้น

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๒๖

๒๗) นายวิรัช โมกแก้ว

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๒๗

๒๘) นายวัชรพงษ์ เทพดนตรี

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๒๘

๒๙) นายอนุศาสน์ สวดี

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๒๙

๓๐) นายกรวิทย์ เจียศิริสกุล

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๓๐

๓๑) นางสาวอริกา รังสวัสดิ์

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๓๑

๓๒) นางสาวนภสวรรณ คงชา

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๓๒

๓๓) นายสุทธิระ อรุณจันทร์

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๓๓

๓๔) นางสาวทัศนีย์ อ่อนคำ

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๓๔

๓๕) นางสาวพริ้มพร

ผู้ยื่นคำขอ...

ผู้อำนวยการกองวิจัยและเตือนภัยมลพิษโรงงาน  
ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม

- ๓๖) นายศุภณัฐ คุณอนกาญจน์  
 ๓๗) นางสาวศิริภาพร เหมือนแร่  
 ๓๘) นางสาวนัส ชานิล  
 ๓๙) นางสาวพรรณิภา ชีระจินดาชล  
 ๔๐) นายนาเคนทร์ พันธุ์ชาติกุล

- ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๓๖  
 ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๓๗  
 ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๓๘  
 ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๓๙  
 ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๔๐



(นางจินดา เดชะศรีพร)

ผู้อำนวยการกองวิจัยและเตือนภัยมลพิษโรงงาน  
 ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม

เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท ยูโนเด็ค แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด เลขทะเบียน ๖-๑๔๕

ที่ อก ๐๓๑๐(๑)/ ๑๘๗ ๙ ลงวันที่ ๐ ๙ กุมภาพันธ์ ๒๕๖๕

ข. เจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๑๐๖ ราย

- |                                |                            |
|--------------------------------|----------------------------|
| ๑) นายสุสันต์ พันสิงห์         | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๐๑ |
| ๒) นางสาวสุธรรมา แก้วช้อนอก    | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๐๒ |
| ๓) นายพริษฐ์ เจริญผล           | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๐๓ |
| ๔) นางสาววิไลลักษณ์ เกโรสง     | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๐๔ |
| ๕) นายสมชาติ อุทุมรัตน์        | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๐๕ |
| ๖) นางสาวปรมาภรณ์ ทองแก้ว      | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๐๖ |
| ๗) นางสาวกัลยา สมพงษ์          | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๐๗ |
| ๘) นายอรรถพร เทพทอง            | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๐๘ |
| ๙) นางสาวอมรรัตน์ พุทธาสี      | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๐๙ |
| ๑๐) นางสาววรรณิ สายบุญเรือน    | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๑๐ |
| ๑๑) นายกฤษณพงษ์ นามทิพย์       | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๑๑ |
| ๑๒) นางสาวอาภรณ์ อ่อนคง        | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๑๒ |
| ๑๓) นายกิตติศักดิ์ ทรงจำรัส    | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๑๓ |
| ๑๔) นางสาวอักษรินทร์ บุญคง     | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๑๔ |
| ๑๕) นางสาวพรพิมล แวนทอง        | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๑๕ |
| ๑๖) นายวิชณุ สุวรรณราช         | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๑๖ |
| ๑๗) นายอภิวิชญ์ ท่วงที         | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๑๗ |
| ๑๘) นายมานิตย์ ปานโชติ         | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๑๘ |
| ๑๙) นายทศพร ธนะพิรุฬห์         | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๑๙ |
| ๒๐) นางสาวกัลยาณี โยธา         | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๒๐ |
| ๒๑) นางสาวเกวลี สุขรี          | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๒๑ |
| ๒๒) นางสาวชมชนัญ อภิพัทธ์ปภา   | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๒๒ |
| ๒๓) นายศิริพัชร จงผดุงเกียรติ  | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๒๓ |
| ๒๔) นางสาวสุภาวดี อินยาศรี     | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๒๔ |
| ๒๕) นายพงศ์เทพ เหล่าขจร        | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๒๕ |
| ๒๖) นายขวัญชัย พันทุกข์        | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๒๖ |
| ๒๗) นางสาวพัชจิรา คดีพิศาล     | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๒๗ |
| ๒๘) นางสาวเมวิกา เสือคำจันทร์  | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๒๘ |
| ๒๙) นายกานต์พงศ์ บุญพวง        | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๒๙ |
| ๓๐) นางสาวพริดา เจริญชัยสมบัติ | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๓๐ |
| ๓๑) นายณพรัตน์ จະโต            | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๓๑ |
| ๓๒) นายพีระพัฒน์ บัญญัติติลป   | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๓๒ |
| ๓๓) นายปรีดา ไชยภูมิสกุล       | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๓๓ |
| ๓๔) นายชัชวาลย์ เลื่อนล่อง     | ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๓๔ |
| ๓๕) นายปิยะณัฐ                 | ๑๔๕-จ-๐๐๓๕                 |

(นางจินดา เดชะศรีพร)

ผู้อำนวยการกองวิจัยและเตือนภัยมลพิษโรงงาน  
 ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม

๓๖) นายณณสินธุ์...

๓๖) นายณกสินธุ์ อนุธรรมรัตน์  
๓๗) นายกันนิก รสโธ  
๓๘) นายจักรพันธ์ ภูมิรินทร์  
๓๙) นายปริญญา กลมเกลียว  
๔๐) นายธีรวัจน์ มาตรโพธิ์ศรี  
๔๑) นายธีรเมธ สุขศรี  
๔๒) นายบุญญฤทธิ์ ก้อนสิน  
๔๓) นายพรธวัช ใวสกุล  
๔๔) นายอชิตะ แสงจันทร์  
๔๕) นายณัฐพงศ์ เมืองชัย  
๔๖) นายธนัท เลิศประเสริฐ  
๔๗) นางสาวนิภาพร จันทเขตต์  
๔๘) นายยุทธพงษ์ อีสระสุข  
๔๙) นายรณภพ ภูตระกูลพัฒนา  
๕๐) นางสาวศิริวรรณ ขอนพา  
๕๑) นายสมพงศ์ สกฤทัย  
๕๒) นายสุริยัน นิธิจิตชูวงศ์  
๕๓) นายอชฎาฐ ยนต์ศิริ  
๕๔) นายเอกภูมิ เสนอใจ  
๕๕) นายสุสันต์ บุญเลี้ยง  
๕๖) นายธนเดช ทวานเสนาะ  
๕๗) นายพิพัฒน์ ต้นธนกุล  
๕๘) นายอภิสิทธิ์ ศรีคงแก้ว  
๕๙) นายภูวดล มงคลสูง  
๖๐) นายอุทัย แก้วรากมุข  
๖๑) นางสาววนารินทร์ สานนท์  
๖๒) นายศุภกร รินวงศ์  
๖๓) นายศักดิ์สิทธิ์ เกิดซัง  
๖๔) นางสาวศิริพร อัมภ์รัตน์  
๖๕) นางสาวจินตสุภา เปี่ยมศรี  
๖๖) นางสาวเนตรนภา กมลบุรณ์  
๖๗) นางสาวอารียา ทารมย์  
๖๘) นายจิรวัฒน์ สุขเกษม  
๖๙) นายกิตติพงษ์ สอนชัยภูมิ  
๗๐) นายจุมพล สอนเพชร  
๗๑) นางสาวพัชราภรณ์ แสงฟ้า  
๗๒) นายรัตนชัย เหล่ามา

ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๓๖  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๓๗  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๓๘  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๓๙  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๐  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๑  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๒  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๓  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๔  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๕  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๖  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๗  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๘  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๙  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๐  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๑  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๒  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๓  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๔  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๕  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๖  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๗  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๘  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๙  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๖๐  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๖๑  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๖๒  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๖๓  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๖๔  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๖๕  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๖๖  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๖๗  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๖๘  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๖๙  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๗๐  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๗๑  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๗๒

๗๓) นายอิทธิพงษ์...

(นางจันทรา เดชะพรหม)  
ผู้อำนวยการกองวิจัยและเตือนภัยมลพิษโรงงาน  
ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม

๗๓) นายอิทธิพงษ์ ศรีวิเศษ  
๗๔) นางสาวกรรณิการ์ สาสีทา  
๗๕) นายธราปกรณ์ พิมพ์ศรี  
๗๖) นายพรชัย คุ้มม่วง  
๗๗) นางสาวทัศนีย์ ไชยหาร  
๗๘) นายธีรพงษ์ ศรีคำแหง  
๗๙) นางสาวณัฐชา พรหมศิริ  
๘๐) นางสาวลัดดาวัลย์ โพธิ์พันธ์  
๘๑) นางสาวกมลวรรณ เจริญจันทร์  
๘๒) นายนพรัตน์ จันทะคุณ  
๘๓) นายปิยวัฒน์ ใหม่ชู  
๘๔) นางสาวพรนัชชา กลิ่นอุณ  
๘๕) นายนภสิทธิ์ ศรีพิมพ์  
๘๖) นางสาวลักขิกา จันทรสุน  
๘๗) นายสงกรานต์ มาลัยทอง  
๘๘) นางสาวสาธิตา แซ่เดียว  
๘๙) นายศักดิ์ศิรินทร์ นุ่มนัม  
๙๐) นายวรพงษ์ นนทจันทร์  
๙๑) นางสาวชนากา มาคะมาตร  
๙๒) นางสาวอนธรัตน์ คุณานุพันธ์ชัย  
๙๓) นายวีระยุทธ สาระภักดิ์  
๙๔) นางสาวธิดิยา วีระพันธุ์วัฒน์  
๙๕) นายกฤตพล พงศ์สถาพร  
๙๖) นายณัฐชัย พรหมอารักษ์  
๙๗) นายชนินทร์ พานแก้ว  
๙๘) นายปรัชชาพล โสภา  
๙๙) นายวัชรินทร์ แสนงาม  
๑๐๐) นางสาวอนภรณ์ ลาพรม  
๑๐๑) นายอาทิตย์ อุดมผล  
๑๐๒) นายปวรร บุนนาค  
๑๐๓) นายอิทธิเดช ใจบุญ  
๑๐๔) นายคณิติน พงษ์อิศรานุพร  
๑๐๕) นางสาวสุภารัตน์ จันทร์ประทีป  
๑๐๖) นายเสฏฐวุฒิ เอ็มกลิ่นบัว

ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๗๓  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๗๔  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๗๕  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๗๖  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๗๗  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๗๘  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๗๙  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๘๐  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๘๑  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๘๒  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๘๓  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๘๔  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๘๕  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๘๖  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๘๗  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๘๘  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๘๙  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๙๐  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๙๑  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๙๒  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๙๓  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๙๔  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๙๕  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๙๖  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๙๗  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๙๘  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๙๙  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๐๐  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๐๑  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๐๒  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๐๓  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๐๔  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๐๕  
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๐๖

(นางจันทรา เดชะพรหม)  
ผู้อำนวยการกองวิจัยและเตือนภัยมลพิษโรงงาน  
ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม

เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท ยูโนเต็ด แอนนาไลซิส แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด เลขทะเบียน ๖-๑๔๕

ที่ อก ๐๓๑๐(๑)/ ๑๘๗๕ ลงวันที่ ๐๙ กุมภาพันธ์ ๒๕๖๕

ขอขยายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๓๕๗ รายการ

น้ำเสีย จำนวน 46 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Aldrin	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
2	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>[4]</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>
3	Barium	Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>
4	$\alpha$ -BHC	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
5	$\beta$ -BHC	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
6	$\delta$ -BHC	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
7	$\gamma$ -BHC	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
8	Biochemical Oxygen Demand	1) 5-Day BOD Test, Azide Modification Method <sup>[4]</sup> 2) 5-Day BOD Test, Membrane Electrode Method <sup>[4]</sup>
9	Cadmium	1) Digestion, Direct Air-Acetylene Flame Method <sup>[4]</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>[4]</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>
10	Chemical Oxygen Demand	1) Closed Reflux, Titrimetric Method <sup>[4]</sup> 2) Closed Reflux, Colorimetric Method <sup>[4]</sup> 3) Open Reflux, Titrimetric Method <sup>[4]</sup>
11	Chlordane	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
12	Chromium	1) Digestion, Direct Air-Acetylene Flame Method <sup>[4]</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>[4]</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>
13	Color	ADMI Weighted-Ordinate Spectrophotometric Method <sup>[4]</sup>
14	Copper	1) Digestion, Direct Air-Acetylene Flame Method <sup>[4]</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>[4]</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>
15	Cyanide	1) Distillation, Colorimetric Method <sup>[4]</sup> 2) Flow Injection Analysis Method <sup>[4]</sup>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
16	o,p'-DDT	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
17	4,4'-DDD	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
18	4,4'-DDE	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
19	4,4'-DDT	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
20	Dieldrin	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
21	Endosulfan I	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
22	Endosulfan II	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
23	Endosulfan sulfate	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
24	Endrin	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
25	Endrin aldehyde	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
26	Formaldehyde	Distillation, Colorimetric Method <sup>[3]</sup>
27	Free Chlorine	1) Iodometric Method <sup>[4]</sup> 2) DPD Ferrous Titrimetric Method <sup>[4]</sup>
28	Heptachlor	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
29	Heptachlor Epoxide	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
30	Hexavalent Chromium	1) Colorimetric Method <sup>[4]</sup> 2) Extraction, Direct Air-Acetylene Flame Method <sup>[4]</sup>
31	Lead	1) Digestion, Direct Air-Acetylene Flame Method <sup>[4]</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>[4]</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>
32	Manganese	1) Digestion, Direct Air-Acetylene Flame Method <sup>[4]</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>[4]</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>
33	Mercury	Digestion, Cold-Vapor Atomic Absorption Spectrometric Method <sup>[4]</sup>
34	Methoxychlor	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
35	Nickel	1) Digestion, Direct Air-Acetylene Flame Method <sup>[4]</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>[4]</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
36	Oil & Grease	1) Liquid-Liquid, Partition-Gravimetric Method <sup>[4]</sup> 2) Soxhlet Extraction Method <sup>[4]</sup>
37	pH	Electrometric Method <sup>[4]</sup>
38	Phenols	1) Distillation, Chloroform Extraction Method <sup>[4]</sup> 2) Distillation, Direct Photometric Method <sup>[4]</sup>
39	Selenium	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>[4]</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>
40	Sulfide	1) Iodometric Method <sup>[4]</sup> 2) Methylene Blue Method <sup>[4]</sup>
41	Temperature	Laboratory and Field Methods <sup>[4]</sup>
42	Total Dissolved Solids	Dried at 180 °C <sup>[4]</sup>
43	Total Kjeldahl Nitrogen	Semi-Micro-Kjeldahl Method <sup>[4]</sup>
44	Total Suspended Solids	Dried at 103-105 °C <sup>[4]</sup>
45	Trivalent Chromium	1) Digestion, Direct Air-Acetylene Flame Method; Colorimetric Method; Calculation <sup>[4]</sup> 2) Digestion, Inductively Coupled Plasma Method; Colorimetric Method; Calculation <sup>[4]</sup>
46	Zinc	1) Digestion, Direct Air-Acetylene Flame Method <sup>[4]</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>[4]</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>

## น้ำใต้ดิน จำนวน 126 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Acenaphthene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
2	Acetone	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
3	Aldrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup> <i>อีก 1 มล</i>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
4	Anthracene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
5	Antimony	Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>
6	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>[4]</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>
7	Atrazine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
8	Barium	1) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>[4]</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>
9	Benz(a)anthracene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
10	Benzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
11	Benzo(b)fluoranthene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
12	Benzo(k)fluoranthene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
13	Benzoic acid	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
14	Benzo(a)pyrene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup> <i>อีก 1 มล</i>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
15	Benzo(g,h,i)perylene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
16	Beryllium	Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>
17	Bis(2-chloroethyl)ether	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
18	Bis(2-ethylhexyl)phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
19	Bromodichloromethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
20	Bromoform	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
21	Butanol	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
22	Butyl benzyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
23	Cadmium	1) Digestion, Direct Air-Acetylene Flame Method <sup>[4]</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>[4]</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>
24	Carbazole	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
25	Carbon disulfide	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
26	Carbon tetrachloride	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
27	Chlordane	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
28	p-Chloroaniline	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
29	Chlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup> <i>อื่นๆ</i>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
30	Chlorodibromomethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
31	Chloroform	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
32	2-Chlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
33	Chromium	1) Digestion, Direct Air-Acetylene Flame Method <sup>[4]</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>[4]</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>
34	Chromium (III)	1) Digestion, Direct Air-Acetylene Flame Method; Colorimetric Method; Calculation <sup>[4]</sup> 2) Digestion, Inductively Coupled Plasma Method; Colorimetric Method; Calculation <sup>[4]</sup>
35	Chromium (VI)	1) Colorimetric Method <sup>[4]</sup> 2) Extraction, Air-Acetylene Flame Method <sup>[4]</sup>
36	Chrysene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
37	Cyanide	Distillation, Colorimetric Method <sup>[4]</sup>
38	2,4-D	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
39	DDD	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
40	DDE	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
41	DDT	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup> <i>อื่นๆ</i>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
42	Dibenz(a,h)anthracene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
43	Di-n-butyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
44	1,2-Dichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
45	1,3-Dichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
46	1,4-Dichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
47	3,3'-Dichlorobenzidine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
48	1,1-Dichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
49	1,2-Dichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
50	1,1-Dichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
51	cis-1,2-Dichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
52	trans-1,2-Dichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
53	2,4-Dichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
54	1,2-Dichloropropane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
55	1,3-Dichloropropane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
56	1,3-Dichloropropene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
57	Dieldrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup> <i>สังเกต</i>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
58	Diethyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
59	2,4-Dimethylphenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
60	2,4-Dinitrophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
61	2,4-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
62	2,6-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
63	Di-n-Octyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
64	Endosulfan	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
65	Endrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
66	Ethylbenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
67	Fluoranthene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
68	Fluorene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
69	Heptachlor	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup> <i>สังเกต</i>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
70	Heptachlor epoxide	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
71	Hexachlorobenzene	2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
72	Hexachloro-1,3-butadiene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
73	n-Hexane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
74	$\alpha$ -HCH	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
75	$\beta$ -HCH	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
76	$\gamma$ -HCH	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
77	Hexachlorocyclopentadiene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
78	Hexachloroethane	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
79	Indeno(1,2,3-cd)pyrene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
80	Isophorone	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
81	Lead	1) Digestion, Direct Air-Acetylene Flame Method <sup>[4]</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>[4]</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
82	Manganese	1) Digestion, Direct Air-Acetylene Flame Method <sup>[4]</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>[4]</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>
83	Mercury	Digestion, Cold-Vapor Atomic Absorption Spectrometric Method <sup>[4]</sup>
84	Methanol	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
85	Methoxychlor	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup>
86	Methyl bromide	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
87	Methylene chloride	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
88	2-Methylphenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
89	2-Methylnaphthalene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
90	Methyl tert-butyl ether	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
91	Naphthalene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
92	Nickel	1) Digestion, Direct Air-Acetylene Flame Method <sup>[4]</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>[4]</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>
93	Nitrobenzene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
94	N-Nitrosodiphenylamine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
95	N-Nitrosodi-n-propylamine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
96	Polychlorinated Biphenyls - PCB 1016 - PCB 1221 - PCB 1232 - PCB-1242 - PCB-1248 - PCB-1254 - PCB-1260	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
97	Pentachlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
98	pH	Electrometric Method <sup>[4]</sup>
99	Phenanthrene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
100	Phenol	1) Distillation, Chloroform Extraction Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
101	Pyrene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
102	Selenium	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>[4]</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>
103	Silver	Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>
104	Styrene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
105	1,1,2,2-Tetrachloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
106	Tetrachloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
107	Toluene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
108	Toxaphene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[4]</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
109	TPH (C <sub>5</sub> - C <sub>8</sub> )	1) Purge and Trap, Gas Chromatographic Method <sup>[11,21]</sup> 2) Purge and Trap, Gas Chromatographic/Mass spectrometric Method <sup>[11,25]</sup>
110	TPH (C <sub>9</sub> - C <sub>16</sub> )	Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[9,21]</sup>
111	TPH (C <sub>16</sub> - C <sub>35</sub> )	Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[9,21]</sup>
112	1,2,4-Trichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
113	1,1,1-Trichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
114	1,1,2-Trichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
115	Trichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
116	2,4,5-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
117	2,4,6-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
118	1,3,5-Trimethylbenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
119	Vanadium	Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>
120	Vinyl acetate	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
121	Vinyl chloride	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
122	m-Xylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
123	o-Xylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
124	p-Xylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
125	Xylene (Total)	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>[4]</sup>
126	Zinc	1) Digestion, Direct Air-Acetylene Flame Method <sup>[4]</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>[4]</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>[4]</sup>

## อากาศเสีย (ปล่อยระบาย) จำนวน 25 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Antimony	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>[5]</sup>
2	Arsenic	1) Isokinetic Sampling, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>[5]</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>[5]</sup>
3	Cadmium	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method <sup>[5]</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>[5]</sup>
4	Carbon Monoxide	Instrumental Analyzer Method <sup>[5]</sup>
5	Chlorine	Isokinetic Sampling, Ion Chromatographic Method <sup>[5]</sup>
6	Chromium	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method <sup>[5]</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>[5]</sup>
7	Cobalt	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>[5]</sup>
8	Copper	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method <sup>[5]</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>[5]</sup>
9	Cresol	Absorption Sampling, Gas Chromatographic Method <sup>[5]</sup>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
10	Dioxins/Furans	Isokinetic Sampling <sup>[5]</sup>
11	Hydrogen Chloride	Isokinetic Sampling, Ion Chromatographic Method <sup>[5]</sup>
12	Hydrogen Fluoride	Isokinetic Sampling, Ion Chromatographic Method <sup>[5]</sup>
13	Hydrogen Sulfide	Absorption Sampling, Iodometric Method <sup>[5]</sup>
14	Lead	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method <sup>[5]</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>[5]</sup>
15	Manganese	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method <sup>[5]</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>[5]</sup>
16	Mercury	Isokinetic Sampling, Digestion, Cold-Vapor Atomic Absorption Spectrometric Method <sup>[5]</sup>
17	Nickel	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method <sup>[5]</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>[5]</sup>
18	Opacity	Ringelmann's Method <sup>[1]</sup>
19	Oxides of Nitrogen	1) Absorption Sampling, Phenoldisulfonic acid Method <sup>[5]</sup> 2) Instrumental Analyzer Method <sup>[5]</sup>
20	Selenium	1) Isokinetic Sampling, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>[5]</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>[5]</sup>
21	Sulfur Dioxide	1) Absorption Sampling, Barium-Thorin Titrimetric Method <sup>[5]</sup> 2) Instrumental Analyzer Method <sup>[5]</sup>
22	Sulfuric Acid	Isokinetic Sampling, Barium-Thorin Titrimetric Method <sup>[5]</sup>
23	Total Suspended Particulate	Isokinetic Sampling, Gravimetric Method <sup>[5]</sup>
24	Vanadium	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>[5]</sup>
25	Xylene	1) Bag Sampling, Gas Chromatographic Method <sup>[5]</sup> 2) Adsorption Sampling, Gas Chromatographic Method <sup>[5]</sup>

## สิ่งปลูกหรือวัสดุที่ไม่ใช่แล้ว จำนวน 35 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Aldrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[2,9,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup>
2	Antimony	Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
3	Arsenic	1) Waste Extraction, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>[2,6,15]</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>[2,6,13]</sup> 3) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>[7,15]</sup> 4) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
4	Barium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>[2,6,13]</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
5	Beryllium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>[2,6,13]</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
6	Cadmium	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method <sup>[2,6,14]</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>[2,6,13]</sup> 3) Digestion, Flame Atomic Absorption Spectrometric Method <sup>[7,14]</sup> 4) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
7	Chlordane	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[2,9,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup>
8	Chromium	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method <sup>[2,6,14]</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>[2,6,13]</sup> <i>เพิ่ม</i>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
9	Chromium (III)	3) Digestion, Flame Atomic Absorption Spectrometric Method <sup>[7,14]</sup> 4) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup> 1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method; Waste Extraction, Colorimetric Method; Calculation <sup>[2,6,14,16]</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method; Waste Extraction, Colorimetric Method; Calculation <sup>[2,6,13,16]</sup> 3) Digestion, Flame Atomic Absorption Spectrometric Method; Alkaline Digestion, Colorimetric Method; Calculation <sup>[7,8,14,16]</sup> 4) Digestion, Inductively Coupled Plasma Method; Alkaline Digestion, Colorimetric Method; Calculation <sup>[7,8,13,16]</sup>
10	Chromium (VI)	1) Waste Extraction, Colorimetric Method <sup>[2,16]</sup> 2) Alkaline Digestion, Colorimetric Method <sup>[8,16]</sup>
11	Cobalt	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>[2,6,13]</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
12	Copper	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method <sup>[2,6,14]</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>[2,6,13]</sup> 3) Digestion, Flame Atomic Absorption Spectrometric Method <sup>[7,14]</sup> 4) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
13	2,4-D	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[2,9,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup>
14	DDD	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[2,9,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup> <i>เพิ่ม</i>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
15	DDE	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[2,9,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup>
16	DDT	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[2,9,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup>
17	Dieldrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[2,9,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup>
18	Endrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[2,9,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup>
19	Heptachlor	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[2,9,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup>
20	Lead	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method <sup>[2,6,14]</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>[2,6,13]</sup> 3) Digestion, Flame Atomic Absorption Spectrometric Method <sup>[7,14]</sup> 4) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
21	Lindane	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[2,9,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup>
22	Mercury	1) Waste Extraction, Digestion, Cold-Vapor Atomic Absorption Spectrometric Method <sup>[2,17]</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>[2,6,13]</sup> <i>อีก</i>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
		3) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method <sup>[18]</sup> 4) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup> 5) Thermal Decomposition Amalgamation and Atomic Absorption Spectrometric Method <sup>[19]</sup>
23	Methoxychlor	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[2,9,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup>
24	Molybdenum	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>[2,6,13]</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
25	Nickel	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method <sup>[2,6,14]</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>[2,6,13]</sup> 3) Digestion, Flame Atomic Absorption Spectrometric Method <sup>[7,14]</sup> 4) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
26	Polychlorinated Biphenyls - Aroclor 1016 - Aroclor 1221 - Aroclor 1232 - Aroclor 1242 - Aroclor 1248 - Aroclor 1254 - Aroclor 1260 - 2-Chlorobiphenyl - 2,3-Dichlorobiphenyl - 2,2',5-Trichlorobiphenyl - 2,4',5-Trichlorobiphenyl - 2,2',3,5'-Tetrachlorobiphenyl - 2,2',5,5'-Tetrachlorobiphenyl - 2,3',4,4'-Tetrachlorobiphenyl - 2,2',3,4,5'-Pentachlorobiphenyl	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[2,9,23]</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,23]</sup> <i>อีก</i>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
27	- 2,2',4,5,5'- Pentachlorobiphenyl	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[2,9,28]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup> Electrometric Method <sup>[31,32]</sup> 3) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>[2,6,20]</sup> 4) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
	- 2,3,3',4',6'- Pentachlorobiphenyl	
	- 2,2',3,4,4',5'- Hexachlorobiphenyl	
	- 2,2',3,4,5,5'- Hexachlorobiphenyl	
	- 2,2',3,5,5',6'- Hexachlorobiphenyl	
	- 2,2',4,4',5,5'- Hexachlorobiphenyl	
	- 2,2',3,3',4,4',5'- Heptachlorobiphenyl	
	- 2,2',3,4,4',5,5'- Heptachlorobiphenyl	
	- 2,2',3,4,4',5,6'- Heptachlorobiphenyl	
	- 2,2',3,4',5,5',6'- Heptachlorobiphenyl	
	- 2,2',3,3',4,4',5,5',6'- Nonachlorobiphenyl	
	Pentachlorophenol	
28	pH	
29	Selenium	

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
30	Silver	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>[2,6,13]</sup>
31	Thallium	2) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup> 1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>[2,6,13]</sup>
32	Toxaphene	2) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup> 1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>[2,9,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup>
33	Trichloroethylene	1) Waste Extraction, Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[2,12,25]</sup> 2) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
34	Vanadium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>[2,6,13]</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
35	Zinc	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method <sup>[2,6,14]</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>[2,6,13]</sup> 3) Digestion, Flame Atomic Absorption Spectrometric Method <sup>[7,14]</sup> 4) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>

## ดิน จำนวน 125 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Acenaphthene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,24]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
2	Acetone	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
3	Aldrin	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
4	Anthracene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,26]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
5	Antimony	Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
6	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>[7,15]</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
7	Atrazine	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
8	Barium	Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
9	Benz(a)anthracene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,24]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
10	Benzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
11	Benzo(b)fluoranthene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,24]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
12	Benzo(k)fluoranthene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,24]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
13	Benzoic acid	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
14	Benzo(a)pyrene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,24]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
15	Benzo(g,h,i)perylene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
16	Beryllium	Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
17	Bis(2-chloroethyl)ether	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
18	Bis(2-ethylhexyl)phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
19	Bromodichloromethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
20	Bromoform	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
21	Butanol	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
22	Butyl benzyl phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
23	Cadmium	1) Digestion, Flame Atomic Absorption Spectrometric Method <sup>[7,14]</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
24	Carbazole	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
25	Carbon disulfide	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
26	Carbon tetrachloride	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
27	Chlordane	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
28	p-Chloroaniline	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
29	Chlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
30	Chlorodibromomethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
31	Chloroform	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
32	2-Chlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
33	Chromium	1) Digestion, Flame Atomic Absorption Spectrometric Method <sup>[7,14]</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
34	Chromium (III)	1) Digestion, Flame Atomic Absorption Spectrometric Method; Alkaline Digestion, Colorimetric Method; Calculation <sup>[7,8,14,16]</sup> 2) Digestion, Inductively Coupled Plasma Method; Alkaline Digestion, Colorimetric Method; Calculation <sup>[7,8,13,16]</sup>
35	Chromium (VI)	Alkaline Digestion, Colorimetric Method <sup>[8,16]</sup>
36	Chrysene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,24]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
37	Cyanide	Extraction, Distillation, Colorimetric Method <sup>[28,29,30]</sup>
38	2,4-D	Ultrasonic Extraction, Gas Chromatographic Method <sup>[27]</sup>
39	DDD	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
40	DDE	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
41	DDT	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
42	Dibenz(a,h)anthracene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,24]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
43	Di-n-butyl phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
44	1,2-Dichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
45	1,3-Dichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
46	1,4-Dichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
47	3,3'-Dichlorobenzidine	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
48	1,1-Dichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
49	1,2-Dichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
50	1,1-Dichloroethylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
51	cis-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
52	trans-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
53	2,4-Dichlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
54	1,2-Dichloropropane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
55	1,3-Dichloropropane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
56	1,3-Dichloropropene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
57	Dieldrin	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
58	Diethyl phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
59	2,4-Dimethylphenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
60	2,4-Dinitrophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
61	2,4-Dinitrotoluene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
62	2,6-Dinitrotoluene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
63	Di-n-Octyl phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
64	Endosulfan	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
65	Endrin	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
66	Ethylbenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
67	Fluoranthene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,24]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
68	Fluorene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,24]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
69	Heptachlor	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
70	Heptachlor epoxide	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
71	Hexachlorobenzene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
72	Hexachloro-1,3-butadiene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
73	n-Hexane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
74	$\alpha$ -HCH	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
75	$\beta$ -HCH	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
76	$\gamma$ -HCH	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
77	Hexachlorocyclopentadiene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
78	Hexachloroethane	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
79	Indeno(1,2,3-cd)pyrene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,24]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
80	Isophorone	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
81	Lead	1) Digestion, Flame Atomic Absorption Spectrometric Method <sup>[7,14]</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
82	Manganese	1) Digestion, Flame Atomic Absorption Spectrometric Method <sup>[7,14]</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
83	Mercury	1) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method <sup>[18]</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup> 3) Thermal Decomposition Amalgamation and Atomic Absorption Spectrometric Method <sup>[19]</sup>
84	Methanol	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
85	Methoxychlor	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
86	Methyl bromide	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
87	Methylene chloride	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
88	2-Methylphenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
89	2-Methylnaphthalene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
90	Methyl tert-butyl ether	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
91	Naphthalene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,24]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
92	Nickel	1) Digestion, Flame Atomic Absorption Spectrometric Method <sup>[7,14]</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
93	Nitrobenzene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
94	N-Nitrosodiphenylamine	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
95	N-Nitrosodi-n-propylamine	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup> 31กพ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
96	Polychlorinated Biphenyls - Aroclor 1016 - Aroclor 1221 - Aroclor 1232 - Aroclor 1242 - Aroclor 1248 - Aroclor 1254 - Aroclor 1260 Polychlorinated Biphenyls - 2-Chlorobiphenyl - 2,3-Dichlorobiphenyl - 2,2',5'-Trichlorobiphenyl - 2,4',5'-Trichlorobiphenyl - 2,2',3,5'-Tetrachlorobiphenyl - 2,2',5,5'-Tetrachlorobiphenyl - 2,3',4,4'-Tetrachlorobiphenyl - 2,2',3,4,5'-Pentachlorobiphenyl - 2,2',4,5,5'-Pentachlorobiphenyl - 2,3,3',4',6-Pentachlorobiphenyl - 2,2',3,4,4',5'-Pentachlorobiphenyl - 2,2',3,4,5,5'-Hexachlorobiphenyl - 2,2',3,5,5',6-Hexachlorobiphenyl - 2,2',4,4',5,5'-Hexachlorobiphenyl - 2,2',3,3',4,4',5-Heptachlorobiphenyl - 2,2',3,4,4',5,5'-Heptachlorobiphenyl - 2,2',3,4,4',5',6-Heptachlorobiphenyl	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,23]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>  Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,23]</sup> 31กพ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
97	- 2,2',3,4',5,5',6-Heptachlorobiphenyl - 2,2',3,3',4,4',5,5',6-Nonachlorobiphenyl Pentachlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
98	Phenanthrene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,24]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
99	Phenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
100	Pyrene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,24]</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
101	Selenium	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>[7,22]</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
102	Silver	Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
103	Styrene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
104	1,1,2,2-Tetrachloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
105	Tetrachloroethylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
106	Toluene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
107	Toxaphene	Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,22]</sup>
108	TPH (C <sub>5</sub> -C <sub>8</sub> )	1) Purge and Trap, Gas Chromatographic Method <sup>[12,21]</sup> 2) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
109	TPH (C <sub>8</sub> -C <sub>16</sub> )	Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,21]</sup>
110	TPH (C <sub>16</sub> -C <sub>35</sub> )	Ultrasonic Extraction, Gas Chromatographic Method <sup>[10,21]</sup>
111	1,2,4-Trichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>

112 1,1,1-Trichloroethane...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
112	1,1,1-Trichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
113	1,1,2-Trichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
114	Trichloroethylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
115	2,4,5-Trichlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
116	2,4,6-Trichlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>[10,26]</sup>
117	1,3,5-Trimethylbenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
118	Vanadium	Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>
119	Vinyl acetate	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
120	Vinyl chloride	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
121	m-Xylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
122	o-Xylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
123	p-Xylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
124	Xylene (Total)	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>[12,25]</sup>
125	Zinc	1) Digestion, Flame Atomic Absorption Spectrometric Method <sup>[7,14]</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>[7,13]</sup>

**เอกสารอ้างอิง**

1. กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2549. เรื่อง กำหนดค่าปริมาณเข้มข้นที่เจือปนในอากาศที่ระบายออกจากปล่องของหม้อน้ำโรงสีข้าวที่ใช้ถ่านเป็นเชื้อเพลิง.

ราชกิจจานุเบกษา. 4 ธันวาคม 2549. เล่มที่ 123 ตอนพิเศษ 125 ง.

2. กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2548. เรื่อง การกำจัดสิ่งปฏิกูลหรือวัสดุที่ไม่ใช้แล้ว. ราชกิจจานุเบกษา. 25 มกราคม 2549. เล่มที่ 123 ตอนพิเศษ 11 ง.

3. สมาคมวิศวกรรมสิ่งแวดล้อมแห่งประเทศไทย. คู่มือวิเคราะห์น้ำเสีย. พิมพ์ครั้งที่ 4. กรุงเทพฯ: เรือนแก้วการพิมพ์, 2547.

4. APHA, AWWA, WEF. **Standard Methods for the Examination of Water and Wastewater**. 23<sup>rd</sup> ed. Washington, DC: APHA, 2017.

5. United States Environmental Protection Agency. **Standards of Performance for New Stationary Sources**. 40 CFR 60. Appendix A, 2019.

6. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods**. SW-846, 1997.

7. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Acid Digestion of Sediments, Sludges, and Soils**. SW-846 Method 3050B, 1996.

8. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Alkaline Digestion for Hexavalent Chromium**. SW-846 Method 3060A, 1996.

9. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste 3. Physical/Chemical Methods. Separatory Funnel Liquid-Liquid Extraction**. SW-846 Method 3510C, 1996.

10. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Ultrasonic Extraction**. SW-846 Method 3550C, 2007.

11. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Purge and Trap for Aqueous Samples**. SW-846 Method 5030C, 2003.

12. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Closed System Purge and Trap and Extraction for Volatile Organics in Soil and Waste Sample**. SW-846 Method 5035A, 2000.

13. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Inductively Coupled Plasma-Optical Emission Spectrometry**. SW-846 Method 6010D, 2014.

14. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Flame Atomic Absorption Spectrophotometry**. SW-846 Method 7000B, 2007.

15. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Arsenic (Atomic Absorption, Gaseous Hydride)**. SW-846 Method 7061A, 1992. *สำเนา*

16. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Chromium, Hexavalent (Colorimetric)**. SW-846 Method 7196A, 1992.

17. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Mercury in Liquid Waste (Manual Cold Vapor Technique)**. SW-846 Method 7470A, 1994.

18. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Mercury in Solid or Semisolid Waste (Manual Cold-Vapor Technique)**. SW-846 Method 7471B, 1998.

19. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Mercury in Solids and Solutions by Thermal Decomposition, Amalgamation, and Atomic Absorption Spectrophotometry**. SW-846 Method 7473, 2007.

20. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Selenium (Atomic Absorption, Borohydride Reduction)**. SW-846 Method 7742, 1994.

21. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Nonhalogenated Organics Using GC/FID**. SW-846 Method 8015D, 2003.

22. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Organochlorine Pesticides by Gas Chromatography**. SW-846 Method 8081B, 2007.

23. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Polychlorinated Biphenyls (PCBs) by Gas Chromatography**. SW-846 Method 8082A, 2007.

24. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Polynuclear Aromatic Hydrocarbons**. SW-846 Method 8100, 1980.

25. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Volatile Organic Compounds by Gas Chromatography/Mass Spectrometry**. SW-846 Method 8260D, 2018.

26. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Semivolatile Organic Compounds by Gas Chromatography/Mass Spectrometry**. SW-846 Method 8270E, 2018.

27. United States Environmental Protection Agency. **Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Chlorinated Herbicides by GC Using Methylation or Pentafluorobenzoylation Derivatization**. SW-846 Method 8151A, 1996. *สำเนา*

28. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. **Total and Amenable Cyanide : Distillation. SW-846 Method 9010C**, 2004.

29. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. **Cyanide Extraction Procedure for Solids and Oils. SW-846 Method 9013A**, 2014.

30. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. **Cyanide in Waters and Extracts using Titrimetric and Manual Spectrophotometric Procedures. SW-846 Method 9014**, 2014.

31. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. **pH Electrometric Measurement. SW-846 Method 9040C**, 2004.

32. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. **Soil and Waste pH. SW-846 Method 9045D**, 2004 *amal*



## ภาคผนวก 14    กฎหมายที่เกี่ยวข้อง

---



## ประกาศกระทรวงทรัพยากรธรรมชาติและสิ่งแวดล้อม

เรื่อง กำหนดมาตรฐานควบคุมการระบายน้ำทิ้ง

จากอาคารบางประเภทและบางขนาด

โดยที่ได้มีการปฏิรูประบบราชการโดยให้มีการจัดตั้งกระทรวงทรัพยากรธรรมชาติและสิ่งแวดล้อมขึ้นมา และให้โอนภารกิจของกระทรวงวิทยาศาสตร์ เทคโนโลยีและสิ่งแวดล้อม ในส่วนที่เกี่ยวข้องกับพระราชบัญญัติส่งเสริมและรักษาคุณภาพสิ่งแวดล้อมแห่งชาติ พ.ศ. ๒๕๓๕ ไปเป็นของกระทรวงทรัพยากรธรรมชาติและสิ่งแวดล้อม ประกอบกับเป็นการสมควรให้คณะกรรมการควบคุมมลพิษเป็นผู้พิจารณาเห็นชอบกับวิธีการตรวจหาค่ามาตรฐานการระบายน้ำทิ้ง นอกเหนือจากวิธีการที่กำหนดไว้ เหนือกรมควบคุมมลพิษ จึงสมควรแก้ไขปรับปรุงประกาศกระทรวงวิทยาศาสตร์ เทคโนโลยีและสิ่งแวดล้อม เรื่อง กำหนดมาตรฐานควบคุมการระบายน้ำทิ้งจากอาคารบางประเภทและบางขนาด

อาศัยอำนาจตามความในมาตรา ๕๕ แห่งพระราชบัญญัติส่งเสริมและรักษาคุณภาพสิ่งแวดล้อมแห่งชาติ พ.ศ. ๒๕๓๕ แก้ไขโดยมาตรา ๑๑๔ แห่งพระราชกฤษฎีกาแก้ไขบทบัญญัติให้สอดคล้องกับการโอนอำนาจหน้าที่ของส่วนราชการ ให้เป็นไปตามพระราชบัญญัติปรับปรุงกระทรวง ทบวง กรม พ.ศ. ๒๕๔๕ พ.ศ. ๒๕๔๕ อันเป็นพระราชบัญญัติที่มีบทบัญญัติบางประการเกี่ยวกับการจำกัดสิทธิและเสรีภาพของบุคคล ซึ่งมาตรา ๒๕ ประกอบกับมาตรา ๓๕ มาตรา ๔๘ มาตรา ๕๐ และมาตรา ๕๑ ของรัฐธรรมนูญแห่งราชอาณาจักรไทยบัญญัติให้กระทำได้ โดยอาศัยอำนาจตามบทบัญญัติแห่งกฎหมาย รัฐมนตรีว่าการกระทรวงทรัพยากรธรรมชาติและสิ่งแวดล้อม โดยคำแนะนำของคณะกรรมการควบคุมมลพิษ และโดยความเห็นชอบของคณะกรรมการสิ่งแวดล้อมแห่งชาติ จึงออกประกาศไว้ ดังต่อไปนี้

ข้อ ๑ ให้ยกเลิกประกาศกระทรวงวิทยาศาสตร์ เทคโนโลยีและสิ่งแวดล้อม เรื่อง กำหนดมาตรฐานควบคุมการระบายน้ำทิ้งจากอาคารบางประเภทและบางขนาด ลงวันที่ ๑๐ มกราคม พ.ศ. ๒๕๓๖

ข้อ ๒ ในประกาศนี้

“อาคาร” หมายความว่า อาคารที่ก่อสร้างขึ้น ไม่ว่าจะมัลักษณะเป็นอาคารหลังเดียว หรือเป็นกลุ่มของอาคารซึ่งตั้งอยู่ในพื้นที่ซึ่งเป็นบริเวณเดียวกัน และไม่ว่าจะมีท่อระบายน้ำท่อเดียว หรือมีหลายท่อที่เชื่อมติดต่อกันระหว่างอาคารหรือไม่ก็ตาม ซึ่งได้แก่

(๑) อาคารชุด ตามกฎหมายว่าด้วยอาคารชุด

(๒) โรงแรม ตามกฎหมายว่าด้วยโรงแรม

(๓) หอพัก ตามกฎหมายว่าด้วยหอพัก

(๔) สถานบริการประเภทสถานอาบน้ำ นวดหรืออบตัว ซึ่งมีผู้ให้บริการแก่ลูกค้า ตามกฎหมายว่าด้วยสถานบริการ

(๕) โรงพยาบาลของทางราชการหรือสถานพยาบาล ตามกฎหมายว่าด้วยสถานพยาบาล

(๖) อาคารโรงเรียนเอกชน ตามกฎหมายว่าด้วยโรงเรียนเอกชน โรงเรียนของทางราชการ อาคารสถานอุดมศึกษาของเอกชน ตามกฎหมายว่าด้วยสถาบันอุดมศึกษาของเอกชนและสถานอุดมศึกษาของทางราชการ

(๗) อาคารที่ทำการของทางราชการ รัฐวิสาหกิจ หรือองค์การระหว่างประเทศและของเอกชน

(๘) อาคารของศูนย์การค้าหรือห้างสรรพสินค้า

(๙) ตลาด ตามกฎหมายว่าด้วยการสาธารณสุข แต่ไม่รวมถึง ท่าเทียบเรือประมง สะพานปลา หรือกิจการแพปลา

(๑๐) ภัตตาคารหรือร้านอาหาร

“น้ำทิ้ง” หมายความว่า น้ำเสียที่ผ่านระบบบำบัดน้ำเสียแล้วจนเป็นไปตามมาตรฐานควบคุมการระบายน้ำทิ้งตามที่กำหนดไว้ในประกาศนี้

ข้อ ๓ ให้แบ่งประเภทของอาคารตามข้อ ๒ ออกเป็น ๕ ประเภท คือ

(๑) อาคารประเภท ก.

(๒) อาคารประเภท ข.

(๓) อาคารประเภท ค.

(๔) อาคารประเภท ง.

(๕) อาคารประเภท จ.

ข้อ ๔ อาคารประเภท ก. หมายความว่า

(๑) อาคารชุดที่มีจำนวนห้องสำหรับใช้เป็นที่อยู่อาศัยรวมกันทุกชั้นของอาคาร หรือกลุ่มของอาคารตั้งแต่ ๕๐๐ ห้องนอนขึ้นไป

(๒) โรงแรมที่มีจำนวนห้องสำหรับใช้เป็นที่พักรวมกันทุกชั้นของอาคาร หรือกลุ่มของอาคารตั้งแต่ ๒๐๐ ห้องขึ้นไป

(๓) โรงพยาบาลของทางราชการ รัฐวิสาหกิจหรือสถานพยาบาล ตามกฎหมายว่าด้วยสถานพยาบาล ที่มีเตียงสำหรับผู้ป่วยไว้ค้างคืนรวมกันทุกชั้นของอาคารหรือกลุ่มของอาคารตั้งแต่ ๓๐ เตียงขึ้นไป



- (๑) หอพักที่มีจำนวนห้องสำหรับใช้เป็นที่อยู่อาศัยรวมกันทุกชั้นของอาคาร หรือกลุ่มของอาคาร ตั้งแต่ ๑๐ ห้อง แต่ไม่ถึง ๕๐ ห้อง
- (๒) ตลาดที่มีพื้นที่ใช้สอยรวมกันทุกชั้นของอาคาร หรือกลุ่มของอาคารตั้งแต่ ๕๐๐ ตารางเมตร แต่ไม่ถึง ๑,๐๐๐ ตารางเมตร
- (๓) ภัตตาคารหรือร้านอาหารที่มีพื้นที่ให้บริการรวมกันทุกชั้นของอาคาร หรือกลุ่มของอาคาร ตั้งแต่ ๑๐๐ ตารางเมตร แต่ไม่ถึง ๒๕๐ ตารางเมตร
- ข้อ ๘ อาคารประเภท จ. หมายความว่า ภัตตาคารหรือร้านอาหารที่มีพื้นที่ให้บริการรวมกันทุกชั้นไม่ถึง ๑๐๐ ตารางเมตร
- ข้อ ๙ มาตรฐานควบคุมการระบายน้ำทิ้งจากอาคาร ประเภท ก. ต้องมีค่าดังต่อไปนี้
- (๑) ความเป็นกรดและด่าง (PH) ต้องมีค่าระหว่าง ๕-๙
- (๒) บีโอดี (BOD) ต้องมีค่าไม่เกิน ๒๐ มิลลิกรัมต่อลิตร
- (๓) สารแขวนลอย (Suspended Solids) ต้องมีค่าไม่เกิน ๓๐ มิลลิกรัมต่อลิตร
- (๔) ซัลไฟด์ (Sulfide) ต้องมีค่าไม่เกิน ๑.๐ มิลลิกรัมต่อลิตร
- (๕) สารที่ละลายได้ทั้งหมด (Total Dissolved Solids) ต้องมีค่าเพิ่มขึ้นจากปริมาณสารละลายในน้ำใช้ตามปกติไม่เกิน ๕๐๐ มิลลิกรัมต่อลิตร
- (๖) ตะกอนหนัก (Settleable Solids) ต้องมีค่าไม่เกิน ๐.๕ มิลลิตรต่อลิตร
- (๗) น้ำมันและไขมัน (Fat Oil and Grease) ต้องมีค่าไม่เกิน ๒๐ มิลลิกรัมต่อลิตร
- (๘) ทีเคเอ็น (TKN) ต้องมีค่าไม่เกิน ๓๕ มิลลิกรัมต่อลิตร
- ข้อ ๑๐ มาตรฐานควบคุมการระบายน้ำทิ้งจากอาคาร ประเภท ข. ต้องเป็นไปตามข้อ ๕ เว้นแต่
- (๑) บีโอดี ต้องมีค่าไม่เกิน ๓๐ มิลลิกรัมต่อลิตร
- (๒) สารแขวนลอย ต้องมีค่าไม่เกิน ๔๐ มิลลิกรัมต่อลิตร
- ข้อ ๑๑ มาตรฐานควบคุมการระบายน้ำทิ้งจากอาคาร ประเภท ค. ต้องเป็นไปตามข้อ ๕ เว้นแต่
- (๑) บีโอดี ต้องมีค่าไม่เกิน ๔๐ มิลลิกรัมต่อลิตร
- (๒) สารแขวนลอย ต้องมีค่าไม่เกิน ๕๐ มิลลิกรัมต่อลิตร

- (๓) ซัลไฟด์ ต้องมีค่าไม่เกิน ๓.๐ มิลลิกรัมต่อลิตร
- (๔) ค่าทีเคเอ็น ต้องมีค่าไม่เกิน ๔๐ มิลลิกรัมต่อลิตร
- ข้อ ๑๒ มาตรฐานควบคุมการระบายน้ำทิ้งจากอาคาร ประเภท ง. ต้องเป็นไปตามข้อ ๕ เว้นแต่
- (๑) บีโอดี ต้องมีค่าไม่เกิน ๕๐ มิลลิกรัมต่อลิตร
- (๒) สารแขวนลอย ต้องมีค่าไม่เกิน ๕๐ มิลลิกรัมต่อลิตร
- (๓) ซัลไฟด์ ต้องมีค่าไม่เกิน ๔.๐ มิลลิกรัมต่อลิตร
- (๔) ค่าทีเคเอ็น ต้องมีค่าไม่เกิน ๔๐ มิลลิกรัมต่อลิตร
- ข้อ ๑๓ มาตรฐานควบคุมการระบายน้ำทิ้งจากอาคาร ประเภท จ. ต้องมีค่าดังต่อไปนี้
- (๑) ความเป็นกรดและด่างต้องมีค่าระหว่าง ๕-๙
- (๒) บีโอดี ต้องมีค่าไม่เกิน ๒๐๐ มิลลิกรัมต่อลิตร
- (๓) สารแขวนลอย ต้องมีค่าไม่เกิน ๖๐ มิลลิกรัมต่อลิตร
- (๔) น้ำมันและไขมัน ต้องมีค่าไม่เกิน ๑๐๐ มิลลิกรัมต่อลิตร
- ข้อ ๑๔ การตรวจสอบมาตรฐานการระบายน้ำทิ้งจากอาคาร ให้ใช้วิธีการดังต่อไปนี้
- (๑) การตรวจสอบค่าความเป็นกรดและด่างให้กระทำโดยใช้เครื่องวัดความเป็นกรดและด่างของน้ำ (PH Meter)
- (๒) การตรวจสอบค่าบีโอดีให้กระทำโดยใช้วิธีการอะไซด์โมดิฟิเคชั่น (Azide Modification) ที่อุณหภูมิ ๒๐ องศาเซลเซียส เป็นเวลา ๕ วัน คิดค่ากันหรือวิธีการอื่นที่คณะกรรมการควบคุมมลพิษให้ความเห็นชอบ
- (๓) การตรวจสอบค่าสารแขวนลอยให้กระทำโดยใช้วิธีการกรองผ่านกระดาษกรองใยแก้ว (Glass Fibre Filter Disc)
- (๔) การตรวจสอบค่าซัลไฟด์ให้กระทำโดยใช้วิธีการไตเตรท (Titrate)
- (๕) การตรวจสอบค่าสารที่ละลายได้ทั้งหมดให้กระทำโดยใช้วิธีการระเหยแห้งระหว่างอุณหภูมิ ๑๐๓ องศาเซลเซียส ถึงอุณหภูมิ ๑๐๕ องศาเซลเซียส ในเวลา ๑ ชั่วโมง
- (๖) การตรวจสอบค่าตะกอนหนักให้กระทำโดยใช้วิธีการกรวยอิมฮอฟฟ์ (Imhoff cone) ขนาดบรรจุ ๑,๐๐๐ ลูกบาศก์เซนติเมตร ในเวลา ๑ ชั่วโมง

(๗) การตรวจสอบค่าน้ำมันและไขมันให้กระทำโดยใช้วิธีการสกัดด้วยตัวทำละลาย แล้วแยกหาน้ำมันของน้ำมันและไขมัน

(๘) การตรวจสอบค่าที่เคเอ็นให้กระทำโดยใช้วิธีการเจลดาล์ (Kjeldahl)

ข้อ ๑๕ การคิดคำนวณพื้นที่ใช้สอย จำนวนอาคารและจำนวนห้องของอาคาร หรือกลุ่มของอาคาร ให้เป็นไปตามวิธีการที่คณะกรรมการควบคุมมลพิษกำหนด โดยประกาศในราชกิจจานุเบกษา

ข้อ ๑๖ วิธีการเก็บตัวอย่างน้ำ ความถี่ และระยะเวลาในการเก็บตัวอย่างน้ำ ให้เป็นไปตามที่คณะกรรมการควบคุมมลพิษกำหนด โดยประกาศในราชกิจจานุเบกษา

ข้อ ๑๗ ประกาศนี้ให้ใช้บังคับตั้งแต่วันถัดจากวันประกาศในราชกิจจานุเบกษาเป็นต้นไป

ประกาศ ณ วันที่ ๗ พฤศจิกายน พ.ศ. ๒๕๔๘  
ยงยุทธ ดิยะไพรัช  
รัฐมนตรีว่าการกระทรวงทรัพยากรธรรมชาติและสิ่งแวดล้อม



## ประกาศคณะกรรมการสิ่งแวดล้อมแห่งชาติ

ฉบับที่ ๘ (พ.ศ. ๒๕๓๗)

ออกตามความในพระราชบัญญัติส่งเสริมและรักษาคุณภาพสิ่งแวดล้อมแห่งชาติ

พ.ศ. ๒๕๓๕

เรื่อง กำหนดมาตรฐานคุณภาพน้ำในแหล่งน้ำผิวดิน

อาศัยอำนาจตามความในมาตรา ๓๒ (๑) แห่งพระราชบัญญัติส่งเสริมและรักษาคุณภาพสิ่งแวดล้อมแห่งชาติ พ.ศ. ๒๕๓๕ คณะกรรมการสิ่งแวดล้อมแห่งชาติประกาศกำหนดมาตรฐานคุณภาพน้ำในแหล่งน้ำผิวดิน ไว้ดังต่อไปนี้

### หมวด ๑

#### บททั่วไป

ข้อ ๑ ในประกาศนี้

“แหล่งน้ำผิวดิน” หมายถึง แม่น้ำ ลำคลอง หนอง บึง ทะเลสาบ อ่างเก็บน้ำ และแหล่งน้ำสาธารณะอื่นๆ ที่อยู่ในพื้นแผ่นดิน ซึ่งหมายความรวมถึงแหล่งน้ำสาธารณะที่อยู่ในพื้นแผ่นดินบนเกาะด้วย แต่ไม่รวมถึงน้ำบาดาล และในกรณีที่แหล่งน้ำนั้นอยู่ติดกับทะเลให้หมายความถึงแหล่งน้ำที่อยู่ในปากแม่น้ำหรือปากทะเลสาบ

ปากแม่น้ำและปากทะเลสาบให้ถือแนวเขตตามที่กรมเจ้าท่ากำหนด

## หมวด ๒

### ประเภทและมาตรฐานคุณภาพน้ำในแหล่งน้ำผิวดิน

ข้อ ๒ ให้แบ่งแหล่งน้ำผิวดินออกเป็น ๕ ประเภทคือ แหล่งน้ำประเภทที่ ๑ แหล่งน้ำประเภทที่ ๒ แหล่งน้ำประเภทที่ ๓ แหล่งน้ำประเภทที่ ๔ และแหล่งน้ำประเภทที่ ๕

(๑) แหล่งน้ำประเภทที่ ๑ ได้แก่ แหล่งน้ำที่คุณภาพน้ำมีสภาพตามธรรมชาติโดยปราศจากน้ำทิ้งจากกิจกรรมทุกประเภทและสามารถเป็นประโยชน์เพื่อ

- (ก) การอุปโภคและบริโภคโดยต้องผ่านการฆ่าเชื้อโรคตามปกติก่อน
- (ข) การขยายพันธุ์ตามธรรมชาติของสิ่งมีชีวิตระดับพื้นฐาน
- (ค) การอนุรักษ์ระบบนิเวศน์ของแหล่งน้ำ

(๒) แหล่งน้ำประเภทที่ ๒ ได้แก่ แหล่งน้ำที่ได้รับน้ำทิ้งจากกิจกรรมบางประเภทและสามารถเป็นประโยชน์เพื่อ

(ก) การอุปโภคและบริโภคโดยต้องผ่านการฆ่าเชื้อโรคตามปกติและผ่านกระบวนการปรับปรุงคุณภาพน้ำทั่วไปก่อน

- (ข) การอนุรักษ์สัตว์น้ำ
- (ค) การประมง
- (ง) การว่ายน้ำและกีฬาทางน้ำ

(๓) แหล่งน้ำประเภทที่ ๓ ได้แก่ แหล่งน้ำที่ได้รับน้ำทิ้งจากกิจกรรมบางประเภทและสามารถเป็นประโยชน์เพื่อ

(ก) การอุปโภคและบริโภคโดยต้องผ่านการฆ่าเชื้อโรคตามปกติและผ่านกระบวนการปรับปรุงคุณภาพน้ำทั่วไปก่อน

- (ข) การเกษตร

(๔) แหล่งน้ำประเภทที่ ๔ ได้แก่ แหล่งน้ำที่ได้รับน้ำทิ้งจากกิจกรรมบางประเภทและสามารถเป็นประโยชน์เพื่อ

(ก) การอุปโภคและบริโภคโดยต้องผ่านการฆ่าเชื้อโรคตามปกติและผ่านกระบวนการปรับปรุงคุณภาพน้ำเป็นพิเศษก่อน

- (ข) การอุตสาหกรรม

(๕) แหล่งน้ำประเภทที่ ๕ ได้แก่ แหล่งน้ำที่ได้รับน้ำทิ้งจากกิจกรรมบางประเภท และสามารถเป็นประโยชน์เพื่อการคมนาคม

ข้อ ๓ คุณภาพน้ำในแหล่งน้ำประเภทที่ ๑ ต้องมีสภาพตามธรรมชาติ และสามารถใช้ประโยชน์ได้ตามข้อ ๒ (๑)

ข้อ ๔ คุณภาพน้ำในแหล่งน้ำประเภทที่ ๒ ต้องมีมาตรฐานดังต่อไปนี้

(๑) ไม่มีวัตถุหรือสิ่งของที่เกิดจากการกระทำของมนุษย์ซึ่งจะทำให้ สี กลิ่น และรสของน้ำเปลี่ยนไปตามธรรมชาติ

(๒) อุณหภูมิ (Temperature) ไม่สูงกว่าอุณหภูมิตามธรรมชาติเกิน ๓ องศาเซลเซียส

(๓) ความเป็นกรดและด่าง (pH) มีค่าระหว่าง ๕.๐-๘.๐

(๔) ออกซิเจนละลาย (DO) มีค่าไม่น้อยกว่า ๖.๐ มิลลิกรัมต่อลิตร

(๕) บีโอดี (BOD) มีค่าไม่เกินกว่า ๑.๕ มิลลิกรัมต่อลิตร

(๖) แบคทีเรียกลุ่มโคลิฟอร์มทั้งหมด (Total Coliform Bacteria) มีค่าไม่เกินกว่า ๕,๐๐๐ เอ็ม.พี.เอ็น. ต่อ ๑๐๐ มิลลิลิตร

(๗) แบคทีเรียกลุ่มฟีคอลโคลิฟอร์ม (Fecal Coliform Bacteria) มีค่าไม่เกินกว่า ๑,๐๐๐ เอ็ม.พี.เอ็น. ต่อ ๑๐๐ มิลลิลิตร

(๘) ไนเตรต ( $\text{NO}_3$ ) ในหน่วยไนโตรเจน มีค่าไม่เกินกว่า ๕.๐ มิลลิกรัมต่อลิตร

(๙) แอมโมเนีย ( $\text{NH}_3$ ) ในหน่วยไนโตรเจน มีค่าไม่เกินกว่า ๐.๕ มิลลิกรัมต่อลิตร

(๑๐) ฟีนอล (Phenols) มีค่าไม่เกินกว่า ๐.๐๐๕ มิลลิกรัมต่อลิตร

(๑๑) ทองแดง (Cu) มีค่าไม่เกินกว่า ๐.๑ มิลลิกรัมต่อลิตร

(๑๒) นิกเกิล (Ni) มีค่าไม่เกินกว่า ๐.๑ มิลลิกรัมต่อลิตร

(๑๓) แมงกานีส (Mn) มีค่าไม่เกินกว่า ๑.๐ มิลลิกรัมต่อลิตร

(๑๔) สังกะสี (Zn) มีค่าไม่เกินกว่า ๑.๐ มิลลิกรัมต่อลิตร

(๑๕) แคดเมียม (Cd) ในน้ำที่มีความกระด้างในรูปของ  $\text{CaCO}_3$  ไม่เกินกว่า ๑๐๐ มิลลิกรัมต่อลิตร มีค่าไม่เกินกว่า ๐.๐๐๕ มิลลิกรัมต่อลิตร และในน้ำที่มีความกระด้างในรูปของ  $\text{CaCO}_3$  เกินกว่า ๑๐๐ มิลลิกรัมต่อลิตร มีค่าไม่เกินกว่า ๐.๐๕ มิลลิกรัมต่อลิตร

(๑๖) โครเมียมชนิดเฮกซะวาเลนต์ (Cr Hexavalent) มีค่าไม่เกินกว่า ๐.๐๕ มิลลิกรัมต่อลิตร

(๑๗) ตะกั่ว (Pb) มีค่าไม่เกิน ๐.๐๕ มิลลิกรัมต่อลิตร

(๑๘)ปรอททั้งหมด (Total Hg) มีค่าไม่เกินกว่า ๐.๐๐๒ มิลลิกรัมต่อลิตร

(๑๙) สารหนู (As) มีค่าไม่เกินกว่า ๐.๐๑ มิลลิกรัมต่อลิตร

(๒๐) ไซยาไนด์ (Cyanide) มีค่าไม่เกินกว่า ๐.๐๐๕ มิลลิกรัมต่อลิตร

(๒๑) กัมมันตภาพรังสี (Radioactivity) มีค่ารังสีแอลฟา (Alpha) ไม่เกินกว่า ๐.๑ เบกเคอเรลต่อลิตร และรังสีเบตา (Beta) ไม่เกินกว่า ๑.๐ เบกเคอเรลต่อลิตร

(๒๒) สารฆ่าศัตรูพืชและสัตว์ชนิดที่มีคลอรีนทั้งหมด (Total Organochlorine Pesticides) มีค่าไม่เกินกว่า ๐.๐๕ มิลลิกรัมต่อลิตร

(๒๓) ดีดีที (DDT) มีค่าไม่เกินกว่า ๑.๐ ไมโครกรัมต่อลิตร

(๒๔) บีเอชซีชนิดแอลฟา (Alpha-BHC) มีค่าไม่เกินกว่า ๐.๐๒ ไมโครกรัมต่อลิตร

(๒๕) ดีลด์ริน (Dieldrin) มีค่าไม่เกินกว่า ๐.๑ ไมโครกรัมต่อลิตร

(๒๖) อัลดริน (Aldrin) มีค่าไม่เกินกว่า ๐.๑ ไมโครกรัมต่อลิตร

(๒๗) เฮปตาคลอร์ (Heptachlor) และเฮปตาคลอร์อีพอกไซด์ (Heptachlorepoxyde) มีค่าไม่เกินกว่า ๐.๒ ไมโครกรัมต่อลิตร

(๒๘) เอนดริน (Endrin) ไม่สามารถตรวจพบได้ตามวิธีการตรวจสอบที่กำหนด

**ข้อ ๕** คุณภาพน้ำในแหล่งน้ำประเภทที่ ๓ ต้องมีมาตรฐานตาม ข้อ ๔ เว้นแต่

(๑) ออกซิเจนละลาย มีค่าไม่น้อยกว่า ๔.๐ มิลลิกรัมต่อลิตร

(๒) บีโอดี มีค่าไม่เกินกว่า ๒.๐ มิลลิกรัมต่อลิตร

(๓) แบคทีเรียกลุ่มโคลิฟอร์มทั้งหมด มีค่าไม่เกินกว่า ๒๐,๐๐๐ เอ็ม.พี.เอ็น.

ต่อ ๑๐๐ มิลลิลิตร

(๔) แบคทีเรียกลุ่มฟีคอลโคลิฟอร์ม มีค่าไม่เกินกว่า ๔,๐๐๐ เอ็ม.พี.เอ็น.

ต่อ ๑๐๐ มิลลิลิตร

**ข้อ ๖** คุณภาพน้ำในแหล่งน้ำประเภทที่ ๔ ต้องมีมาตรฐานตามข้อ ๔ (๑) ถึง (๕) และ (๘) ถึง (๒๘) เว้นแต่

(๑) ออกซิเจนละลาย มีค่าไม่น้อยกว่า ๒.๐ มิลลิกรัมต่อลิตร

(๒) บีโอดี มีค่าไม่เกินกว่า ๔.๐ มิลลิกรัมต่อลิตร

**ข้อ ๗** คุณภาพน้ำในแหล่งน้ำประเภทที่ ๕ ต้องมีมาตรฐานต่ำกว่าคุณภาพน้ำ ในแหล่งน้ำประเภทที่ ๔

**ข้อ ๘** การกำหนดให้แหล่งน้ำผิวดินแหล่งใดแหล่งหนึ่งเป็นประเภทใดตามข้อ ๒ ให้เป็นไปตามที่กรมควบคุมมลพิษประกาศในราชกิจจานุเบกษา

### หมวด ๓

#### วิธีการเก็บตัวอย่างและตรวจสอบคุณภาพน้ำในแหล่งน้ำผิวดิน

**ข้อ ๙** การเก็บตัวอย่างน้ำเพื่อตรวจสอบคุณภาพตามข้อ ๓ ถึง ข้อ ๗ ให้ใช้วิธีการดังต่อไปนี้

(๑) แหล่งน้ำไหล ซึ่งได้แก่ แม่น้ำ ลำคลอง เป็นต้น ให้เก็บที่จุดกึ่งกลางความกว้างของแหล่งน้ำที่ระดับกึ่งกลางความลึก ณ จุดตรวจสอบ เว้นแต่แบคทีเรียกลุ่มโคลิฟอร์มทั้งหมดและแบคทีเรียกลุ่มฟีคอลโคลิฟอร์ม ให้เก็บที่ระดับความลึก ๓๐ เซนติเมตร ณ จุดตรวจสอบ

(๒) แหล่งน้ำนิ่ง ซึ่งได้แก่ ทะเลสาบ หนอง บึง อ่างเก็บน้ำ เป็นต้น ให้เก็บที่ระดับความลึก ๑ เมตร ณ จุดตรวจสอบสำหรับแหล่งน้ำที่มีความลึกเกินกว่า ๒ เมตร และให้เก็บที่จุดกึ่งกลางความลึก ณ จุดตรวจสอบสำหรับแหล่งน้ำที่มีความลึกไม่เกิน ๒ เมตร เว้นแต่แบคทีเรียกลุ่มโคลิฟอร์มทั้งหมดและแบคทีเรียกลุ่มฟีคอลโคลิฟอร์ม ให้เก็บที่ระดับความลึก ๓๐ เซนติเมตร ณ จุดตรวจสอบ

จุดตรวจสอบตาม (๑) และ (๒) ของแหล่งน้ำที่กำหนดตามข้อ ๘ ให้เป็นไปตามที่กรมควบคุมมลพิษกำหนด

**ข้อ ๑๐** การตรวจสอบคุณภาพน้ำตามข้อ ๓ ถึงข้อ ๗ ให้ใช้วิธีการดังต่อไปนี้

(๑) การตรวจสอบอุณหภูมิ ให้ใช้เครื่องวัดอุณหภูมิ (Thermometer) วัดขณะทำการเก็บตัวอย่างน้ำ

(๒) การตรวจสอบค่าความเป็นกรดและด่าง ให้ใช้เครื่องวัดความเป็นกรดและด่างของน้ำ (pH meter) ตามวิธีการหาค่าแบบอิเล็กโตรเมตริก (Electrometric)

(๓) การตรวจสอบค่าออกซิเจนละลาย ให้ใช้วิธีอะไซด์โมดิฟิเคชัน (Azide Modification)

(๔) การตรวจสอบค่าบีไอดี ให้ใช้วิธีอะไซด์โมดิฟิเคชัน (Azide Modification) ที่อุณหภูมิ ๒๐ องศาเซลเซียส เป็นเวลา ๕ วันติดต่อกัน

(๕) การตรวจสอบค่าแบคทีเรียกลุ่มโคลิฟอร์มทั้งหมดและค่าแบคทีเรียกลุ่มฟิคอลโคลิฟอร์ม ให้ใช้วิธีมัลติเพิล ทิวบ์ เฟอว์เมนเตชัน เทคนิค (Multiple Tube Fermentation Technique)

(๖) การตรวจสอบค่าไนเตรตในหน่วยไนโตรเจน ให้ใช้วิธีแคดเมียมรีดักชัน (Cadmium Reduction)

(๗) การตรวจสอบค่าแอมโมเนียในหน่วยไนโตรเจน ให้ใช้วิธีดิสทิลเลชันเนสสเลอร์ไรเซชัน (Distillation Nesslerization)

(๘) การตรวจสอบค่าฟีนอล ให้ใช้วิธีดิสทิลเลชัน ๔ - อะมิโนแอนติไพรีน (Distillation, 4-Amino antipyrine)

(๙) การตรวจสอบค่าทองแดง นิกเกิล แมงกานีส สังกะสี แคดเมียมโครเมียมชนิดเฮกซะวาเลนต์ และตะกั่ว ให้ใช้วิธีอะตอมมิก แอ็บซอร์ปชัน ไดเรกต์ แอสไพเรชัน (Atomic Absorption - Direct Aspiration)

(๑๐) การตรวจสอบค่าปรอททั้งหมด ให้ใช้วิธีอะตอมมิก แอ็บซอร์ปชัน โคลด์เวปเปอร์ เทคนิค (Atomic Absorption-Cold Vapour Technique)

(๑๑) การตรวจสอบค่าสารหนู ให้ใช้วิธีอะตอมมิก แอ็บซอร์ปชัน แก๊สไฮไดรด์ (Atomic Absorption - Gaseous Hydride)

(๑๒) การตรวจสอบค่าไพริไดน์ ให้ใช้วิธีไพริดีน บาร์บิบูริก แอซิด (Pyridine - Barbituric Acid)

(๑๓) การตรวจสอบค่ากัมมันตภาพรังสี ให้ใช้วิธีโลว์ แบ็กกราวด์พร็อพอร์ชันนอล เคาน์เตอร์ (Low Background Proportional Counter)

(๑๔) การตรวจค่าสารฆ่าศัตรูพืชและสัตว์ชนิดที่มีคลอรีนทั้งหมด ดีดีที บีเอชซีชนิดแอลฟา ดีลดริน อัลดริน เฮปตาคลอโรอีพอกไซด์ และเอนดริน ให้ใช้วิธีแก๊สโครมาโตกราฟี (Gas - Chromatography)

**ข้อ ๑๑** การตรวจสอบค่าออกซิเจนละลายให้ใช้ค่าเปอร์เซ็นต์ไทล์ที่ ๒๐ (20<sup>th</sup> Percentile Value) ส่วนการตรวจสอบค่าบีไอดี แบคทีเรียกลุ่มโคลิฟอร์มทั้งหมด และแบคทีเรียกลุ่มฟิคอลโคลิฟอร์ม ให้ใช้ค่าเปอร์เซ็นต์ไทล์ที่ ๘๐ โดยจำนวนและระยะเวลาสำหรับการเก็บตัวอย่างน้ำดังกล่าว ให้เป็นไปตามที่กรมควบคุมมลพิษกำหนด

**ข้อ ๑๒** การเก็บตัวอย่างน้ำตามข้อ ๘ และการตรวจสอบคุณภาพน้ำตามข้อ ๑๐ จะต้องเป็นไปตามวิธีการมาตรฐานสำหรับการวิเคราะห์น้ำและน้ำเสีย (Standard Methods for Examination of Water and Wastewater) ซึ่ง American Public Health Association และ American Water Works Association กับ Water Pollution Control Federation ของสหรัฐอเมริกา ร่วมกันกำหนดไว้ด้วย

ประกาศ ณ วันที่ ๒๐ มกราคม พ.ศ. ๒๕๓๗

ชวน หลีกภัย

นายกรัฐมนตรี

ประธานคณะกรรมการสิ่งแวดล้อมแห่งชาติ

(ประกาศในราชกิจจานุเบกษา เล่ม ๑๑๑ ตอนที่ ๑๖ ง วันที่ ๒๔ กุมภาพันธ์ ๒๕๓๗)