

ภาคผนวก ง

ใบแสดงการตรวจเทียบเครื่องมือ



## High Volume TSP & PM-10 Calibration Data Sheet

Calibration Location : SECOT Co.,Ltd. Calibration Date : Jan 14, 2022  
 Hi-Vol Pump No. : BH-008 Indicator No. : CM-01  
 Amb. Temp (°C) : 25 Press (mmHg) : 760  
 Calibration by : Mr.Punkawin K.

Plate	Indicate (X) ( cm. )	True H <sub>2</sub> O ( in. )	Actual Flow (Y) (cfm)	XY	X <sup>2</sup>	Remark
18	17.40	12.50	58.84	1,023.82	302.76	
13	14.40	10.10	53.20	766.08	207.36	
10	11.60	7.80	46.90	544.04	134.56	
7	7.60	5.10	38.17	290.09	57.76	
5	4.80	3.10	30.04	144.19	23.04	
Sum	55.80	38.60	227.15	2,768.22	725.48	

Calibrated by : Ranbari Approved by : Mr.Haya K.



## High Volume TSP & PM-10 Calibration Data Sheet

Calibration Location : SECOT Co.,Ltd. Calibration Date : Jan 13, 2022  
 Hi-Vol Pump No. : BH-018 Indicator No. : CM-01  
 Amb. Temp (°C) : 25 Press (mmHg) : 760  
 Calibration by : Mr.Punkawin K.

Plate	Indicate (X) ( cm. )	True H <sub>2</sub> O ( in. )	Actual Flow (Y) (cfm)	XY	X <sup>2</sup>	Remark
18	17.00	12.40	58.61	996.37	289.00	
13	14.00	10.20	53.45	748.30	196.00	
10	11.00	8.00	47.48	522.28	121.00	
7	7.20	5.20	38.53	277.42	51.84	
5	4.20	3.20	30.50	128.10	17.64	
Sum	53.40	39.00	228.57	2,672.47	675.48	

Calibrated by : Ranbari Approved by : Mr.Haya K.



## High Volume TSP & PM-10 Calibration Data Sheet

Calibration Location : SECOT Co.,Ltd. Calibration Date : Feb 3, 2022  
 Hi-Vol Pump No. : BH-033 Indicator No. : CM-01  
 Amb. Temp (°C) : 25 Press (mmHg) : 760  
 Calibration by : Mr.Punkawin K.

Plate	Indicate (X) ( cm. )	True H <sub>2</sub> O ( in. )	Actual Flow (Y) (cfm)	XY	X <sup>2</sup>	Remark
18	16.80	12.60	59.07	945.20	256.00	
13	14.00	10.20	53.45	748.30	196.00	
10	11.70	7.90	47.19	528.60	125.40	
7	7.20	5.10	38.17	274.80	51.80	
5	4.00	3.10	30.04	120.20	16.00	
Sum	53.70	38.90	227.92	2,617.10	645.20	

Calibrated by : Punkawin Approved by : Witaya K.



## High Volume TSP & PM-10 Calibration Data Sheet

Calibration Location : SECOT Co.,Ltd. Calibration Date : Feb 3, 2022  
 Hi-Vol Pump No. : BH-035 Indicator No. : CM-01  
 Amb. Temp (°C) : 25 Press (mmHg) : 760  
 Calibration by : Mr.Punkawin K.

Plate	Indicate (X) ( cm. )	True H <sub>2</sub> O ( in. )	Actual Flow (Y) (cfm)	XY	X <sup>2</sup>	Remark
18	19.20	12.80	59.53	1,142.98	368.64	
13	15.40	10.10	53.20	819.28	237.16	
10	12.40	8.10	47.77	592.35	153.76	
7	8.20	5.20	38.53	315.95	67.24	
5	4.80	3.10	30.04	144.19	23.04	
Sum	60.00	39.30	229.07	3,014.74	849.84	

Calibrated by : Punkawin Approved by : Witaya K.

## CERTIFICATE OF ANALYSIS

### Grade of Product: EPA Protocol

Part Number: E04N199E15AC084 Reference Number: 82-401409170-1  
Cylinder Number: EB0102326 Cylinder Volume: 144.4 CF  
Laboratory: 124 - Riverton (SAP) - NJ Cylinder Pressure: 2015 PSIG  
PGVP Number: B52019 Valve Outlet: 660  
Gas Code: CO,NO,NOX,SO2,BALN Certification Date: Feb 05, 2019

Expiration Date: Feb 05, 2027

Certification performed in accordance with "EPA Traceability Protocol for Assay and Certification of Gaseous Calibration Standards (May 2012)" document EPA 600/R-12/531, using the assay procedures listed. Analytical Methodology does not require correction for analytical interference. This cylinder has a total analytical uncertainty as stated below with a confidence level of 95%. There are no significant impurities which affect the use of this calibration mixture. All concentrations are on a volume/volume basis unless otherwise noted.

Do Not Use This Cylinder below 100 psig, i.e. 0.7 megapascals.

ANALYTICAL RESULTS					
Component	Requested Concentration	Actual Concentration	Protocol Method	Total Relative Uncertainty	Assay Dates
NOX	50.00 PPM	51.01 PPM	G1	+/- 0.9% NIST Traceable	01/28/2019, 02/05/2019
NITRIC OXIDE	50.00 PPM	50.86 PPM	G1	+/- 0.9% NIST Traceable	01/28/2019, 02/05/2019
SULFUR DIOXIDE	50.00 PPM	50.87 PPM	G1	+/- 1.0% NIST Traceable	01/28/2019, 02/05/2019
CARBON MONOXIDE	0.5000 %	0.5050 %	G1	+/- 0.7% NIST Traceable	01/31/2019
NITROGEN	Balance				

CALIBRATION STANDARDS					
Type	Lot ID	Cylinder No	Concentration	Uncertainty	Expiration Date
NTRM	13060206	CC401947	4950 PPM CARBON MONOXIDE/NITROGEN	+/- 0.4%	Feb 15, 2019
PRM	12367	APEX1099237	9.82 PPM NITROGEN DIOXIDE/AIR	+/- 2.0%	Jun 02, 2017
NTRM	12010724	KAL004497	50.03 PPM NITRIC OXIDE/NITROGEN	+/- 0.8%	Mar 12, 2024
GMIS	1114201601	CC506710	4.871 PPM NITROGEN DIOXIDE/NITROGEN	+/- 2.0%	Nov 14, 2019
NTRM	14010327	KAL004376	48.08 PPM SULFUR DIOXIDE/NITROGEN	+/- 1.0%	Apr 17, 2024

The SRM, PRM or RGM noted above is only in reference to the GMIS used in the assay and not part of the analysis.

ANALYTICAL EQUIPMENT		
Instrument/Make/Model	Analytical Principle	Last Multipoint Calibration
Siemens Ultramat 6 J3-599 COHIGH	NDIR	Jan 18, 2019
Nicolet 6700 APW1100391 NO	FTIR	Jan 10, 2019
Nicolet 6700 APW1100391 NO2	FTIR	Jan 10, 2019
Nicolet 6700 APW1100391 SO2	FTIR	Jan 10, 2019

Triad Data Available Upon Request

PERMANENT NOTES: PRODUCED IN ACCORDANCE WITH ISO17025 REQUIREMENTS

#### NOTES:

Gross Weight: 27806.3 grams

Net Weight: 4733.2 grams

This calibration std. has been certified in accordance with the May 2012 EPA Traceability Protocol Document EPA-600/R-12/531. All testing processes and measurements conform to the requirements of ISO/IEC 17025 and to Airgas ISO 9001:2008 and relate only to items identified on this certificate. All items are certified to be NIST Traceable with total uncertainty as detailed under Analytical Uncertainty. This document shall not be reproduced in full without written approval of the issuer.



TESTING CERT No. 3082.05

Approved for Release

Page 1 of 82-401409170-1

Sheet No. : CAL-PI-PS10-01/2022



## PITOT TUBE CALIBRATION

Calibration Location: SECOT

Calibration Date: 14/01/2022

Calibrated duct No.: 1

Calibration Standard Pitot tube data

Pitot No.: Sid-01

Coefficient (Cp): 1

Type S Pitot No.: PS10-01

Calibrated by: Mr. Montri P.

### A Side Calibration

Run No.	$\Delta P_{std}$ (mm H <sub>2</sub> O)	$\Delta P_s$ (mm H <sub>2</sub> O)	Cp(s)	Deviation, $\delta$ Cp(s) - Cp(A)
1	7.55	10.75	0.8380	0.0032
2	7.55	10.75	0.8380	0.0032
3	7.55	11.00	0.8285	-0.0064

$C_{P(A),avg}$  0.8349

### B Side Calibration

Run No.	$\Delta P_{std}$ (mm H <sub>2</sub> O)	$\Delta P_s$ (mm H <sub>2</sub> O)	Cp(s)	Deviation, $\delta$ Cp(s) - Cp(B)
1	7.55	11.00	0.8285	-0.0097
2	7.55	10.75	0.8380	-0.0001
3	7.55	10.50	0.8480	0.0098

$C_{P(B),avg}$  0.8382

$|Cp(A) - Cp(B)| = 0.0033$

$C_{P(Avg)} = 0.8365$

Approved by:   
(Miss Katesarin Vorradetwittaya)

\*\*\*  $\delta$  must be  $\leq 0.01$  for the test to be acceptable \*\*\*  
\*\*\*  $|Cp(A) - Cp(B)|$  must also be  $< 0.01$  if average of  $Cp(A)$  and  $Cp(B)$  is to be used \*\*\*



## CONTROL UNIT CALIBRATION

(Metric units, mm)

Date 25 Jan 22

Initial Final Average  
Barometric press, Pb 758 758 758 mmHg

## Dry Gas Meter Data

Console No. M50-06

Metering System ID

DGM Number 333249

DGM Model ES-110

Calibrated by : Montri P.

## Reference Dry Gas Meter Data

Serial No. 358794

Model S110

Correction factor (Yr) 0.9966

Last Calibration Date 8 Jan 22

Orifice manometer setting, ΔH mm H2O	Ref. DGM Volume V <sub>r</sub> Liters	DGM Volume V <sub>m</sub> Liters	Temperature (°C)				Time Θ min	DGM Correction factor (Y)	ΔH@ mm
			Ref DGM T <sub>r</sub>	Dry Gas Meter					
				Inlet T <sub>i</sub>	Outlet T <sub>o</sub>	Avg T <sub>m</sub>			
12.5	100.1	100.6	24	24	23	23.5	8.58	0.9887	42.5446
25.0	100.2	100.2	24	24	23	23.5	6.00	0.9921	41.5532
50.0	100.1	99.7	24	24	23	23.5	4.32	0.9941	43.1019
76.0	100.1	100.9	24	24	23	23.5	3.52	0.9805	43.4295
100.0	100.2	99.6	24	24	23	23.5	3.52	0.9904	42.9584
150.0	100.2	100.5	24	24	23	23.5	2.47	0.9784	42.0708

Average 0.9874 42.6097

Approved by :   
( Miss Katesarin Vorradetwittaya )



## ELECTRICAL AND ELECTRONICS INSTITUTE FOUNDATION FOR INDUSTRIAL DEVELOPMENT

975 Moo 4, Bangpoo Industrial Estate, Soi 8, Sukhumvit Road km 37,

Phraek Sa, Mueang Samut Prakan, Samut Prakan 10280

Tel: +66 2709 4860-8 Fax: +66 2324 0917-8



Certificate No.: CP20210098EA

Operation No.: CP2021120019

## Certificate of Calibration

Equipment: Sound Calibrator

Manufacturer: Cirrus Research Plc

Model/Type: CR:515

Serial No.: 94296

ID No.:

Customer: SECOT Co.,Ltd.


Address: 239 Rimklongprapa Rd., Bangsue,  
Bangkok 10800 Thailand

Received Date: 21 December 2021

Calibrated Date: 24 December 2021

Issued Date: 28 December 2021

Calibrated by: Ms. Juntaporn Kunhakom

Approved by:   
( Mr. Sittichai Swaksuriyawong )  
Group Manager

This report was prepared electronically using applicable electronic signature. Printing or copy of file are considered as a copy of the document.

The reported uncertainty of measurement was based on standard uncertainty multiplied by a coverage factor  $k = 2.00$ , providing a level of confidence of approximately 95%. This certificate may not be reproduced other than in full except with the prior written approval of the Electrical and Electronics Institute, Foundation for Industrial Development.



ELECTRICAL AND ELECTRONICS INSTITUTE  
FOUNDATION FOR INDUSTRIAL DEVELOPMENT

Certificate No.: CP20210098EA

Calibration Report

Equipment: Sound Calibrator  
Manufacturer: Cirrus Research Plc  
Model/Type: CR515  
Serial No.: 94296  
ID No.:  
Ambient Temperature: ( 23 ± 2 ) °C  
Relative Humidity: ( 50 ± 15 ) %  
Pressure: (101.3 ± 1.5) kPa

Method of Calibration :-

IEC 60942:2017

Condition of this result of calibration

1. Reference standards instrument :-

Instrument	Model	Serial No.	Cert. No.	Due Date
1) Standard microphone	4180	2661000	AA-1010-21	13 June 2022
2) Waveform Generator	335118	MY52302264	0144RF21	17 June 2022
3) Audio Analyzing DMM	2015-P	4079144	E1U210398	2 February 2022
4) Pressure humidity and Temperature Transmitter	PTU301	F0640002	CL1-P210047 0255TE21	16 June 2022 7 July 2022

2. This result of calibration was found accurate as shown on date and place of calibration only.

3. This certification is traceable to the international system of unit maintained at :-

Reference standards instrument for Acoustic function

- National Institute of Metrology (Thailand)

Reference standards instrument for Electrical function

- Electrical and Electronics Institute; ONSC Accredited Calibration No.0119

Result of Calibration:-

1. Function : Sound pressure level

Norminal Frequency (Hz)	Specified Sound Pressure level (dB)	Measured value (dB)	Deviated value <sup>[1]</sup> (dB)	Acceptance limit <sup>[3]</sup> (dB)
1000	94	93.80	-0.20	±0.25

2. Function : Frequency

Norminal Sound Pressure level (dB)	Specified Frequency (Hz)	Measured value (Hz)	Deviated value <sup>[2]</sup> (%)	Acceptance limit <sup>[3]</sup> (%)
94	1000	1000.3	0.0	±0.7



ELECTRICAL AND ELECTRONICS INSTITUTE  
FOUNDATION FOR INDUSTRIAL DEVELOPMENT

Certificate No.: CP20210098EA

Calibration Report

3. Function : Total distortion + noise

Norminal Sound Pressure level (dB)	Norminal Frequency (Hz)	Measured value <sup>[4]</sup> (%)	Acceptance limit <sup>[5]</sup> (%)
94	1000	1.4	2.5

Uncertainty of measurement

Function	Uncertainty	Maximum-permitted uncertainty of measurement
Sound pressure level	0.10 dB	0.15 dB
Frequency	0.10 %	0.20 %
Total distortion + noise	0.40 %	0.50 %

Note: [1] The deviated value is the absolute value of the difference between the measured value and the corresponding specified sound pressure level.

[2] The deviated value is the absolute value of the difference in percent between the measured value and the corresponding specified frequency.

[3] The acceptance limit is for the deviated value.

[4] The measured value is the total distortion + noise, measured over the frequency range from 20 Hz to 20 kHz.

[5] The acceptance limit is for the Measured value.

Remarks: 1. Acceptance limit was IEC 60942:2017 Class 1.

- - End of Report - -

# CERTIFICATE OF CALIBRATION

ISSUED BY Noisemeters  
DATE OF ISSUE 06/04/22 CERTIFICATE NUMBER 172691

NoiseMeters

NoiseMeters  
Acoustic House  
Bridlington Road  
Hunmanby  
YO14 0PH  
United Kingdom  
www.noisemeters.com

Page 1 of 1

Test engineer:  
Nigel Smith  
Electronically signed:

*Nigel Smith*

## doseBadge Reader

### Instrument

Manufacturer: Cirrus Research plc  
Model Number: RC:110A  
Serial Number: 95167  
Notes:

### Calibration Procedure

The tests were carried out in accordance with the requirements of IEC 60942:2003 where applicable.

Date of Calibration: 06 April 2022

### Functionality Results

Function	Result
Keypad	Pass
Battery Power	Pass
Display	Pass
Communication	Pass
2 way IR link	Pass
Clock	Pass

### Calibration Results

	Level (dB)	Frequency (Hz)	Distortion (% THD + Noise)
Initial	113.86	990.4	0.48
Adjusted	114.00	990.4	0.48
Uncertainty	± 0.11	± 0.14	± 0.10
Tolerances	± 0.60	± 2.00	± 4.00

### Environmental Conditions

Pressure: 98.30 kPa  
Temperature: 22.8 °C  
Humidity: 40.3 %

### Notes

This certificate provides traceability of measurement to the SI system of units and/or to units of measurement realised at the National Physical Laboratory or other recognised national metrology institutes. This certificate may not be reproduced other than in full, except with the prior written approval of the issuing laboratory. The results within this certificate relate only to the items calibrated. The reported expanded uncertainty is based on a standard uncertainty multiplied by a coverage factor k=2, providing a coverage probability of approximately 95%.



Request Service No. 098/65

Page 1 of 3

## Calibration Certificate

Nomenclature: Brand : Mettler Toledo Type : Top-Loading Electronic Balance

Model : AG245 Serial No. : 1117293916 (198129-0)

Submitted by : Laboratory of SECOT CO., LTD.

Location of Calibration : BAL Room , 6<sup>th</sup> Floor, Secot Co., Ltd.

Calibration range : 0 - 200 g Scale division : 0.00001 g (41g) / 0.0001 g (210g)

Calibration date : May 26, 2022

Reference Standard No. M2110188S, M210183, M220177

Traceable to : Metrological Center SCI ECO Services Co., Ltd, THAI CALIBRATION SERVICES Co., Ltd

Ambient Condition : Temperature 24.28 - 24.42 °C

Humidity 48.10 - 50.90 % RH

Calibrated By : *Sasipa Jaidee* Approved By : *Siripa Jhannong*

( Miss Sasipa Jaidee )

( Miss Siripa Jhannong )

Testing Officer

Chief of Technical Management

Date : 26/05/2022

Date : 26/05/2022

Issued Date : May 27, 2022

## Measurement Report

Request Service No. 098/65

Page 2 of 3

Description: Brand : Mettler Toledo Type : Top-Loading Electronic Balance

Model : AG245                      Serial No. : 1117293916 (198129-0)

Calibration range : 0 - 200 g                      Scale division : 0.00001 g (41g) / 0.0001 g (210g)

Calibration date : May 26,2022

Ambient Condition : Temperature 24.28–24.42 °C Relative humidity 48.10–50.90 % RH

Measurement data :

### 1. Repeatability of Reading :

Load (g)	Standard Deviation of Reading (g)	Maximum Difference between Successive Reading (g)
50	0.000047	0.0002
100	0.000067	0.0002
150	0.000048	0.0001
200	0.000052	0.0001

### 2. Off-Center Loading :

A Mass of 50.0000 g was placed and moved to various position on the pan.

Unit : g

Center	Front	Left	Back	Right	Center	Maximum Difference
50.00020	50.00046	50.00030	50.00000	50.00010	50.00020	0.00026

Issued Date : May 27,2022

Request Service No.098/65

Page 3 of 3

### 3. Departure from Nominal Valve

Reading (g)	Correction (g)	Uncertainty (+/- g)
0	0.00000	$\pm 0.00008$
0.5	0.00005	$\pm 0.00014$
1	-0.00014	$\pm 0.00018$
10	-0.00071	$\pm 0.00034$
20	-0.00091	$\pm 0.00047$
40	-0.00151	$\pm 0.00074$
60	-0.00030	$\pm 0.00012$
80	-0.00021	$\pm 0.00014$
100	-0.00038	$\pm 0.00016$
120	-0.00041	$\pm 0.00018$
140	-0.00048	$\pm 0.00021$
160	-0.00050	$\pm 0.00023$
180	-0.00067	$\pm 0.00025$
200	-0.00124	$\pm 0.00027$

Calibrated by : ..... Susipa ..... Indee .....

(Miss Sasipa Jaidee)

Testing Officer

Date : 26/05/2022

.. Approved By

(Miss Siripa Jhannong)

Chief of Technical Management

Date: 26/05/2022

Issued Date : May 27, 2022





Request Service No.100/65

Page 1 of 3

### Calibration Certificate

Nomenclature : Brand : Sartorius Type : Top-Loading Electronic Balance

Model : BSA224S-CW Serial No. : 32191636

Submitted by : Laboratory of SECOT CO., LTD.

Location of Calibration : BAL Room , 6<sup>th</sup> Floor, Secot Co., Ltd.

Calibration range : 0 – 200 g Scale division : 0.0001 g (220 g)

Calibration date : May 24, 2022

Reference Standard No. M220177, M210183

Traceable to : Metrological Center SCI ECO Services Co., Ltd.

Ambient Condition : Temperature 24.80-24.90 °C

Humidity 50.4-52.9 % RH

Calibrated By : *Khemchuda Insora* Approved By : *Miss Siripa Jhannong*

(Miss Khemchuda Insora)

(Miss Siripa Jhannong)

Testing Officer

Chief of Technical Management

Date : *25/05/2022*

Date : *25/05/2022*

Issued Date : May 25, 2022

### Measurement Report

Request Service No.100/65

Page 2 of 3

Description : Brand : Sartorius

Type : Top-Loading Electronic Balance

Model : BSA224S-CW

Serial No. : 32191636

Calibration range : 0 – 200 g

Scale division : 0.0001 g (220 g)

Calibration date : May 25, 2021

Ambient Condition : Temperature 24.80-24.90 °C Relative humidity 50.4-52.9 % RH

Measurement data :

1. Repeatability of Reading :

Load (g)	Standard Deviation of Reading (g)	Maximum Difference between Successive Reading (g)
50	0.00010	0.0003
100	0.00008	0.0003
150	0.00005	0.0001
200	0.00005	0.0001

2. Off-Center Loading :

A Mass of 50.0000 g was placed and moved to various position on the pan.

Unit : g

Center	Front	Left	Back	Right	Center	Maximum Difference
49.99980	49.99984	49.99994	49.99986	49.99994	49.99980	0.00014

Issued Date : May 25, 2022

## 3. Departure from Nominal Value :

Reading (g)	Correction (g)	Uncertainty (+/- g)
0	0.00000	$\pm 0.00008$
1	-0.00004	$\pm 0.00008$
5	+0.00013	$\pm 0.00008$
10	+0.00018	$\pm 0.00008$
20	+0.00009	$\pm 0.00010$
40	-0.00005	$\pm 0.00010$
60	+0.00012	$\pm 0.00014$
80	+0.00017	$\pm 0.00014$
100	-0.00020	$\pm 0.00017$
120	+0.00003	$\pm 0.00019$
140	+0.00004	$\pm 0.00021$
160	+0.00006	$\pm 0.00022$
180	+0.00004	$\pm 0.00025$
200	+0.00002	$\pm 0.00027$

Calibrated by : Khemchuda Insorn

(Miss Khemchuda Insorn)

Testing Officer

Date : 15/05/2022Approved By : [Signature]

(Miss Siripa Jhannong)

Chief of Technical Management

Date : 25/05/2022

Issued Date : May 25, 2022



## Calibration Certificate

Nomenclature : Brand : Mettler Toledo Type : Top-Loading Electronic Balance

Model : AB204-S Serial No. : 1123163292 (209359)

Submitted by : Laboratory of SECOT CO., LTD.

Location of Calibration : BAL Room , 6<sup>th</sup> Floor, Secot Co., Ltd.

Calibration range : 0 – 200 g Scale division : 0.0001 g (220 g)

Calibration date : May 26, 2022

Reference Standard No. M220177, M210183

Traceable to : Metrological Center SCI ECO Services Company Limited.

Ambient Condition : Temperature 24.80 - 25.30 °C

Humidity 54.6 - 57.5 % RH

Calibrated By : Janista Kui-on

(Miss Janista Kui-on)

Testing Officer

Date : 27/05/2022Approved By : [Signature]

(Miss Siripa Jhannong)

Chief of Technical Management

Date : 27/05/2022

Issued Date : May 27, 2022

## Measurement Report

Request Service No. 101/65

Page 2 of 3

Description : Brand : Mettler Toledo      Type : Top-Loading Electronic Balance  
Model : AB204-S      Serial No. : 1123163292 (209359)  
Calibration range : 0 – 200 g      Scale division : 0.0001 g (220 g)  
Calibration date : May 26, 2022  
Ambient Condition : Temperature 24.80-25.30 °C      Relative humidity 54.6-57.5 % RH

### Measurement data :

#### 1. Repeatability of Reading :

Load (g)	Standard Deviation of Reading (g)	Maximum Difference between Successive Reading (g)
50	0.00008	0.0002
100	0.00010	0.0003
150	0.00010	0.0003
200	0.00012	0.0004

#### 2. Off-Center Loading :

A Mass of 50.0000 g was placed and moved to various position on the pan.

Unit : g

Center	Front	Left	Back	Right	Center	Maximum Difference
50.00004	49.99998	50.00010	50.00008	50.00002	50.00006	0.00006

Issued Date : May 27, 2022

Request Service No.101/65

Page 3 of 3

#### 3. Departure from Nominal Value :

Reading (g)	Correction (g)	Uncertainty (+/- g)
0	0.00000	± 0.00010
1	+0.00004	± 0.00010
5	+0.00011	± 0.00010
10	+0.00002	± 0.00010
20	-0.00005	± 0.00014
40	-0.00010	± 0.00020
60	-0.00004	± 0.00014
80	-0.00012	± 0.00014
100	-0.00024	± 0.00016
120	-0.00025	± 0.00019
140	-0.00020	± 0.00022
160	-0.00022	± 0.00025
180	-0.00015	± 0.00025
200	-0.00016	± 0.00029

Calibrated by : Janista Kui-on

Approved By : 

(Miss Janista Kui-on)

(Miss Siripa Jhannong)

Testing Officer

Chief of Technical Management

Date : 27/05/2022

Date : 27/05/2022

Issued Date : May 27, 2022



## Certificate of Calibration

### ICS-1000 : Cation & Anion (ID#057)

This certificate is to verify that instrument below are calibrated

by Archemica Lab Co.,Ltd.

ICS-1000 S/N : 04090295

RFC-30 S/N : 04100666

for

SECOT Co.,Ltd.



Operator Signature : \_\_\_\_\_ Date : Dec 16, 2021

(Mr.Channarong Khiao-un)

Test Engineer



## BECTHAI BANGKOK EQUIPMENT & CHEMICAL CO., LTD. CALIBRATION LABORATORY

300 Phaholyothin Road, Phayathai, Bangkok 10400, Thailand Tel: +66 2615-2929 Fax: +66 2615-2350-1  
E-mail: bkk@becthai.com Website: www.becthai.com



Certificate No. : CAL-22-119

Page : 1 of 4

## CERTIFICATE OF CALIBRATION

Equipment	:	Spectrophotometer
Manufacturer	:	Thermo Scientific
Model	:	Genesys 150 UV-VIS
Serial No.	:	9A5Y332022
ID No.	:	N/A
Customer	:	Secot Company Limited
	:	239 Rimklongprapa Road,
	:	Bangsue, Bangkok 10800, Thailand
Location	:	Customer Laboratory
Date of Receipt	:	1 March 2022
Date of Calibration	:	1 March 2022
Date of Issue	:	5 March 2022
Ambient Temperature	:	(25±10) °C
Relative Humidity	:	(60±20) %
Condition As-Received	:	Used Item

Calibrated by

( Mr. Palawat Lunchak )

Calibration Engineer

Approved by

( Ms. Jintana Sangthajaroenlap )

Calibration Manager

The reported expanded uncertainty of measurement was based on a combined standard uncertainty multiplied by a coverage factor  $k=2.00$ , providing a level of confidence of approximately 95%.

This certificate may not be reproduced other than in full, except with the prior written approval of the head of Calibration Laboratory.

Indicated values are valid for the state of the Spectrophotometer at the time of calibration only.

ISSUE: 5 REV:4

FM-CAL-33/2

15/05/61



**BECTHAI BANGKOK EQUIPMENT & CHEMICAL CO., LTD.**  
**CALIBRATION LABORATORY**

300 Phaholyothin Road, Phayathai, Bangkok 10400, Thailand Tel: +66 2615-2929 Fax: +66 2615-2350-1  
 E-mail: bkk@becthai.com Website: www.becthai.com



Certificate No. : CAL-22-119

Page : 2 of 4

## CALIBRATION REPORT

### Conditions of this result of calibration

#### 1. Reference Standard Material :

Material	Model	Serial No.	Cert.No.	Due date
Holmium Glass Filter	RM-HG	12705	98236	12 Feb 24
Didymium Glass Filter	RM-DG	13498	98233	12 Feb 24
Neutral Density Filter	RM-1N2N3N	8323	98259	13 Feb 24
Potassium Dichromate Solution	RM-06	23429	98252	12 Feb 24

2. Traceability : This certification is traceable to the International System of Unit maintained at;  
 The Starna Scientific Ltd. Accredited Calibration Laboratory No. 0659.

#### 3. Method of calibration :

The calibration procedure was carried out according to the Guide to CPM-CAL-02 based on ASTM E275-08 (2013) and-  
 ASTM E925-09 (2014).

#### 4. Result of calibration :

( ✓ ) without adjustment

( ) after adjustment

#### 5. Equipment Specifications:

Spectral Bandwidth :	2.0	nm
Data Interval :	0.2	nm
Scan Speed :	Slow	nm/min

FM-CAL-33/2

*Signature* 15/05/61



**BECTHAI BANGKOK EQUIPMENT & CHEMICAL CO., LTD.**  
**CALIBRATION LABORATORY**

300 Phaholyothin Road, Phayathai, Bangkok 10400, Thailand Tel: +66 2615-2929 Fax: +66 2615-2350-1  
 E-mail: bkk@becthai.com Website: www.becthai.com



Certificate No. : CAL-22-119

Page : 3 of 4

## CALIBRATION REPORT

### Wavelength Calibration

Certified Values of Reference Material (nm)	Nominal Value (nm)	UUC*Reading (nm)	Error (nm)	Uncertainty of Measurement (± nm)
241.74	241.74	241.260	-0.480	0.16
637.98	637.98	637.547	-0.433	0.17
879.27	879.27	878.895	-0.375	0.17

### Photometric Calibration for Visible

Wavelength (nm)	Certified Values of Reference Material (A)	UUC* Reading (A)	Error (A)	Uncertainty of Measurement (± A)
420.0	Zero	0.000	0.0000	0.0028
	0.5716	0.571	-0.0006	0.0031
	0.7358	0.733	-0.0028	0.0034
	1.0713	1.072	0.0007	0.0032
440.0	Zero	0.000	0.0000	0.0028
	0.561	0.561	0.0000	0.0031
	0.718	0.716	-0.0020	0.0031
	1.0459	1.047	0.0011	0.0032
465.0	Zero	0.000	0.0000	0.0028
	0.5111	0.511	-0.0001	0.0032
	0.6618	0.660	-0.0018	0.0030
	0.9635	0.965	0.0015	0.0034
546.1	Zero	0.000	0.0000	0.0028
	0.5222	0.521	-0.0012	0.0034
	0.6687	0.667	-0.0017	0.0029
	0.9768	0.977	0.0002	0.0043
590.0	Zero	0.000	0.0000	0.0028
	0.5541	0.553	-0.0011	0.0033
	0.6975	0.695	-0.0025	0.0031
	1.0206	1.020	-0.0006	0.0045
635.0	Zero	0.000	0.0000	0.0028
	0.5398	0.539	-0.0008	0.0033
	0.6658	0.664	-0.0018	0.0033
	0.9741	0.974	-0.0001	0.0045

Remark : Each individual filter is measured against the empty filter holder (blank) used to zero the Spectrophotometer.

Note:

UUC\* : Unit Under Calibration

ISSUE: 5 REV:4

FM-CAL-33/2

*Signature* 15/05/61

ISSUE: 5 REV:4



**BECTHAI BANGKOK EQUIPMENT & CHEMICAL CO., LTD.**  
**CALIBRATION LABORATORY**

300 Phaholyothin Road, Phayathai, Bangkok 10400, Thailand Tel: +66 2615-2929 Fax: +66 2615-2350-1  
 E-mail: bkk@becthai.com Website: www.becthai.com



Certificate No. : CAL-22-119

Page : 4 of 4

## CALIBRATION REPORT

Photometric Calibration for UV

Wavelength (nm)	Certified Values of Reference Material (A)	UUC* Reading (A)	Error (A)	Uncertainty of Measurement ( $\pm$ A)
235.0	Zero	0.000	0.0000	0.0050
	0.7345	0.735	0.0005	0.0076
257.0	Zero	0.000	0.0000	0.0050
	0.8498	0.849	-0.0008	0.0074
313.0	Zero	0.000	0.0000	0.0050
	0.2853	0.287	0.0017	0.0055
350.0	Zero	0.000	0.0000	0.0050
	0.6306	0.630	-0.0006	0.0063

Remark : The Potassium Dichromate Filled cells are measured against a Perchloric acid blank.

Note:

UUC\* : Unit Under Calibration

- End of Report -

ISSUE: 5 REV:4

FM-CAL-33/2

*Signature* 15/05/61



ศูนย์บริการสอบเทียบเพื่ออุตสาหกรรม  
 ศูนย์บริการห้องปฏิบัติการอุตสาหกรรม  
 Foundation for Industrial Development: National Food Institute  
 Food Industrial Laboratory Service Center



## Calibration Certificate

**Certificate No.:** 2203876-001-01  
**Client name:** SECOT CO., LTD.  
**Address:** 239 Rimklongprapa Road,  
 Bangsue, Bangsue, Bangkok 10800

Page 1 of 3

**Equipment:** CHAMBER (Hot Air Oven)  
**Manufacturer:** BINDER  
**Model:** ED 53  
**Serial No.:** 01-27152  
**ID No.:** N/A  
**Order No.:** 2203876  
**Operation No.:** 2203876-001  
**Date of Receipt:** 1 August 2022  
**Date of Calibration:** 1 August 2022

**Calibrated by** Mr.Yothin Charoensuk  
 Scientist

**Approved by** *Signature*  
 (Mr.Pheraphat Tuanjit) (for)  
 Manager, Division of Calibration Laboratory

**Date of Issue:** 3 August 2022

Responsible for the Technical Management Team

The uncertainties are for a confidence probability of approximately 95 %.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65

2008 อาคารอุตสาหกรรม 35 ถนนพหลโยธิน แขวงสามยุค เขตปทุมธานี 10700  
 2008 So. 35, Anon Amarin Road, Bang Yai Khin Subdistrict, Bang Phai District, Bangkok 10700, Thailand  
 Tel: +66(0) 2422 5533 Fax: +66(0) 2422 5545



## Calibration Report

**Certificate No.:** 2203876-001-01  
**Equipment:** CHAMBER (Hot Air Oven)  
 Model: ED 53 Serial No.: 01-27152  
 Resolution: 1 °C ID No.: N/A  
 Manufacturer: BINDER  
**Date of Calibration:** 1 August 2022

Page 2 of 3

**Location:** Laboratory, SECOT CO., LTD.  
**Environment Condition:** Ambient Temperature ( 30 ± 1 ) °C  
 Relative Humidity ( 66 ± 5 ) %  
 Line Voltage ( 220 ± 5 ) Volt

### Condition of this results of Calibration:

- This instrument was calibrated by insert 9 standard thermometer into its chamber and calibration according to W-TE-014 Based on TLAS G-20-1/02-08 (E): Guidelines for Calibration and Checks of Temperature Controlled Enclosures.  
 - The temperature scale used was based on ITS - 90.  
 - All data show below were final values and the initial data may be obtained upon request.

### 2. Reference Standard Instrument :


Instrument	Model	Serial No./ID No.	Certificate No.	Due Date	Through
Digital Thermometer with sensor	34972A	MY57003188	TE 650469-01	11 June 2023	NATIONAL FOOD INSTITUTE
	RTD	CH#101-109/ RTD#101-109			

- This certificate is traceable to International System of Units (SI Units).
- This certificate was certified only for the instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.
- Condition of Calibrated item : Good

### UUC Description :

Time of Record 1 Hour 9 Minute At 104,110 and 180 °C  
 Fresh air Damper ☐ Open Position ☐  
☒ Close Fan ☐  
☐ Not Available

7. Result of Calibration : ☒ Without adjustment ☐ After adjustment

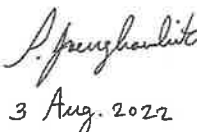  
 3 Aug. 2022



## Calibration Report

**Certificate No.:** 2203876-001-01  
**Equipment:** CHAMBER (Hot Air Oven)  
 Model: ED 53 Serial No.: 01-27152  
 Resolution: 1 °C ID No.: N/A  
 Manufacturer: BINDER

**Date of Calibration:** 1 August 2022

Page 3 of 3

**Calibration point:** 104,110 and 180 °C

### Calibration result:

Calibration Condition	Temperature (°C)	Relative Humidity (%)	Line Voltage (Volt)
MIN	29.8	61	215.0
MAX	30.9	71	225.0

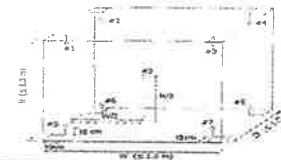


Table 1 : Reporting of Temperature

Calibration point (°C)	Measured Temperature (°C) @ Sensor No. (Sensor No.9 is REF)									Uncertainty ± (°C)
	# 1	# 2	# 3	# 4	# 5	# 6	# 7	# 8	# 9	
104	103.88	104.38	104.57	104.17	103.06	102.86	103.29	103.14	102.94	0.80
110	109.86	110.37	110.58	110.15	109.05	108.83	109.31	109.16	108.93	0.81
180	179.86	180.90	180.31	180.22	179.43	179.49	179.88	180.20	179.67	0.90

Table 2 : Reporting of Characterization Result

UUC* Setting (°C)	UUC* reading (°C)			Stability ± (°C)	Uniformity (°C)	Overall Variation (°C)
	MIN	MAX	Average			
104	104	104	104	0.17	1.6	2.0
110	110	110	110	0.21	1.7	2.0
177	177	177	177	0.33	1.2	2.2

**Note** The quoted uncertainty include " Stability " and " Loading effect (20% of Temp Uniformity) "

UUC\* = Unit Under Calibration


Stability = One-half of the greatest maximum difference of measured temperatures at any one sensors, for at least half an hour after reaching steady state.

Uniformity = The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time.

Overall Variation = The difference of the maximum and minimum measured temperatures throughout observation time.

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k= 2, providing a level of confidence of approximately 95 %.

----- End -----

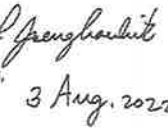  
 3 Aug. 2022



# CAL

Calibratech Co.,Ltd.

7/106-7 Moo 2, Sukhprachasan 3 Rd., Bangpood, Pakkred, Nonthaburi 11120

Tel.(02) 964-6211 Fax.(02) 964-5155, e-mail : calibratech.cal@yahoo.com, calibratech.cal@hotmail.com



NSG-TIS-TIS17025  
CALIBRATION 0030

## Certificate of Calibration

Certificate No. : 65-420016-1

Page : 1 of 2

Submitted by : Secot Co.,Ltd.

239 RimKlongprapa Road, Bangsue, Bangkok 10800 Thailand

Equipment : pH Meter with electrode

pH meter

Manufacturer : Mettler Toledo Model : Seven2Go S2

Range : N/A pH Resolution : 0.01 pH

Serial No. : B924795409 ID No. : N/A

Electrode

Model : InLab Expert Go-ISM Serial No. : 7861180

Environment : Ambient Temperature :  $(25 \pm 2) ^\circ\text{C}$

Relative Humidity :  $(50 \pm 15) \%$

Date of Received : 15 February 2022

Date of Calibration : 24 February 2022

Date of Issue : 24 February 2022

Calibrated by : Bunjerd Masri

Calibration Method : In-house method CAL-M4201 direct measurement by using standard voltage calibrator and using certified reference material (CRM)

Reference Standard Instruments : This certification is traceable to the International System of Units

1. Multiproduct Calibrator

ID No. Cert. No. Due Date Traceability

440001 21E997 17 Mar 2023 National Institute of Metrology Thailand (NIMT)

2. Standard Buffer Solution

pH Cert. No. Lot No. Exp. Date Traceability

4.008 61235182 795894 14 Feb 2024 CPA Chem Ltd. Accredited to ISO 17034 and ISO/IEC 17025

6.985 61223875 769927 15 May 2022 CPA Chem Ltd. Accredited to ISO 17034 and ISO/IEC 17025

10.008 61244986 795895 25 Feb 2023 CPA Chem Ltd. Accredited to ISO 17034 and ISO/IEC 17025

Approved by :

(Bunjerd Masri)

Supervisor

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full except with the prior written approval of the Calibratech Co.,Ltd.



# CAL

Calibratech Co.,Ltd.

7/106-7 Moo 2, Sukhprachasan 3 Rd., Bangpood, Pakkred, Nonthaburi 11120

Tel.(02) 964-6211 Fax.(02) 964-5155, e-mail : calibratech.cal@yahoo.com, calibratech.cal@hotmail.com

## Certificate of Calibration

Certificate No. : 65-420016-1

Page : 2 of 2

Result of Calibration :

UUC Condition As-Received : Good

Function : Electrical measurement

pH meter

Performing standard curve by Multiproduct Calibrator at pH (4,7,10)

Adjustment Curve at nominal pH	Applied Voltage ( mV )	Nominal Value ( pH )	UUC Reading		Correction ( mV )	Uncertainty ( $\pm$ mV )
			( pH )	( mV )		
4, 7, 10	177.4800	4	4.00	177	0	0.58
	0.0000	7	7.00	0	0	0.58
	-177.4800	10	10.00	-177	0	0.58

Function : pH meter with electrode

Performing a three - buffer standard curve using buffer nominal pH (4,7,10)

Adjustment Curve at nominal pH	Standard Buffer ( pH )	UUC Reading ( pH )	Correction ( pH )	Uncertainty ( $\pm$ pH )
4, 7, 10	4.008	4.01	0.00	0.010
	6.985	7.00	-0.01	0.011
	10.008	10.01	0.00	0.014

Remark

UUC : Unit Under Calibration

This result of calibration was found accurate as shown on date and place of calibration only.

This reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor  $k = 2$ , providing a level of confidence of approximately 95%

- 000 -

Bunjerd Masri



CAL-F0031-03





TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES  
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250  
TEL. 0-2717-3000-27 FAX. 0-2719-9484



Cert.No.: 22CH1624  
Page.: 1 of 2

## Certificate of Calibration

Equipment : Conductivity Meter  
Manufacturer : Hanna  
Model : HI98192  
Serial No. : 05200045101  
ID No. : -  
Condition As-Received: Used Item  
Received Date : 22 November 2022  
Calibration Date : 23 November 2022  
Reference : 2211-0761DN-2  
Submitted by : Secot Co., Ltd.  
239 Rimkiongprapa Road,  
Bangsue, Bangkok 10800  
Ambient Temperature : (25  $\pm$  2.5) °C  
Relative Humidity : (50  $\pm$  15) %  
Calibration Procedure: In-house method :  
- CP-CH6 : based on direct measurement by  
using certified reference material (CRM)

Calibrated by : Walalak Sirithean

Approved by :

*Malee*  
Approved Signatory

(✓) Malee Butkruea  
( ) Saithip Meangmai  
( ) Warakorn Lerngagtrakul

Issue Date : 25 November 2022

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Cert.No.: 22CH1624

Page.: 2 of 2

### Condition of this result of calibration

#### 1. Reference Standard Instrument :-

Instrument	Serial No.	ID No.	Certificate No.	Due date
1) Thermometer	9549224	130RC003	221484	17 Apr 2023

This certification is traceable to the International System of Unit maintained at:-

- Traceable to National Institute of Metrology (Thailand), NIMT

#### 2. Certified Reference Materials :-

- Conductivity calibration solution; CPA chem Ltd., The measurement results are traceable to SI through CPA chem Ltd., ANSI-ASQ National Accreditation Board, Accredited No. AR-1835  
- Conductivity calibration solution, Thermo Scientific (traceable to NIST)

Conductivity Solution	Manufacturer	Lot No.	Exp. date
*100 $\mu$ S/cm	Thermo Scientific	152/01	14 Apr 2023
1.413 mS/cm	CPA Chem	823328	20 June 2023
12.880 mS/cm	CPA Chem	823329	20 June 2023

- Control Conductivity calibration solution temperature by Water bath (25 $\pm$ 0.1) °C

3. This certificate is valid only to the item calibrated on date and place of calibration.

### Calibration results

Function : Conductivity Measurement

(\*) After Adjustment at 1.413, 12.880 mS/cm

Conductivity Electrode Serial No.: 0720001N

Standard Conductivity Solution	Before Adjustment UUC* Reading	After Adjustment UUC* Reading	Uncertainty of Measurement ( $\pm$ )	Coverage factor k
*100 $\mu$ S/cm	99.36 $\mu$ S/cm	106.6 $\mu$ S/cm	5.1 $\mu$ S/cm	2.00
1.413 mS/cm	1.296 mS/cm	1.412 mS/cm	0.0093 mS/cm	2.00
12.880 mS/cm	10.86 mS/cm	12.88 mS/cm	0.086 mS/cm	2.00

Remark - UUC\* = Unit Under Calibration  
- \* = Not NSC - ONSC Accredited

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-o0o-

*Malee*  
a 1137186

A 0047740



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)

CALIBRATION AND TESTING EQUIPMENT SERVICES

534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250

TEL. 0-2717-3000-24 FAX. 0-2719-9484

Certificate of Calibration

Equipment : Turbidity Meter  
Manufacturer : Hanna  
Model : HI 98703  
Serial No. : 03090023991  
ID. No. :  
Condition As-Received: Used Item  
Received Date : 6 May 2022  
Calibration Date : 10 May 2022  
Reference : 2205-0123DN-1  
Submitted by : Secot Co., Ltd.  
239 Rimklongprapa Road,  
Bangsue, Bangkok 10800  
Ambient Temperature : (25 ± 2.5) °C  
Relative Humidity : (50 ± 20) %  
Calibration Procedure : In - house method : CP-CH11  
based on direct measurement by  
using Formazin standard solution

Calibrated by : Saithip Meangmai

Approved by :

Approved Signatory

- (✓) Malee Butkruea  
( ) Saithip Meangmai  
( ) Warakorn Lerngagtrakul

Issue Date : 12 May 2022

The Uncertainties are for a confidence probability of approximately 95%.

This certificate may not be reproduced other than in full, except with the prior written approval of the head of Calibration and Testing Equipment Services.

Cert.No.: 22CH633

Page.: 1 of 2



Cert.No. : 22CH633

Page. : 2 of 2

Condition of this calibration result

1. Reference Standard Instruments :

This certification is traceable to the International System of unit (SI unit) through Technology Promotion Association (Thailand-Japan).

Instruments	Serial No.	ID No.	Certificate No.	Due date
1) Thermo-Hygrograph	1103328	130EC010	21H1462	27 June 2022
2) Electronic Balance	N03679	140RC001	21MM429	21 Sep 2022

2. Standard Material : The Formazin suspension has been prepared gravimetric from

Material	Manufacturer	Lot No.	Assay
1) Hexamethylenetetramine	HIMEDIA	0000493947	99.65%
2) Hydrazinium Sulfate	HIMEDIA	0000522014	99.40%

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration result

Performing four - Formazin suspension standard curve by using 0, 20, 100, 800 NTU  
Turbidity Meter Serial Number : 03090023991

Standard Formazine suspension ( NTU )	UUC* Reading ( NTU )	Uncertainty of Measurement ( ± NTU )	Coverage Factor k
20	20.2	0.490	2.05
40	40.0	0.39	2.00
100	100	0.79	2.00
400	400	1.5	2.11

Remark

- UUC\* = Unit Under Calibration  
- NTU = Nephelometric Turbidity Units

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor  $k$ , providing a level of confidence of approximately 95 %.

-o0o-

A 0008849

a 1107001


## Calibration Certificate

**Certificate No.:** 2203876-002-01  
**Client name:** SECOT CO., LTD.  
**Address:** 239 Rimklongprapa Road,  
Bangsue, Bangsue, Bangkok 10800

Page 1 of 3

**Equipment:** CHAMBER (Incubator)  
**Manufacturer:** MEMMERT  
**Model:** ICP 400  
**Serial No.:** K406.0004  
**ID No.:** N/A  
**Order No.:** 2203876  
**Operation No.:** 2203876-002  
**Date of Receipt:** 1 August 2022  
**Date of Calibration:** 1 August 2022

**Calibrated by** Mr.Yothin Charoensuk  
Scientist

**Approved by**   
(Mr.Pheraphat Tuanjit) (for)  
Manager, Division of Calibration Laboratory  
Responsible for the Technical Management Team

**Date of Issue:** 3 August 2022

The uncertainties are for a confidence probability of approximately 95 %.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65



## Calibration Report

**Certificate No.:** 2203876-002-01  
**Equipment:** CHAMBER (Incubator)  
Model: ICP 400 Serial No.: K406.0004  
Resolution: 0.1 °C ID No.: N/A  
Manufacturer: MEMMERT  
**Date of Calibration:** 1 August 2022

Page 2 of 3

**Location:** Laboratory, SECOT CO., LTD.  
**Environment Condition:** Ambient Temperature ( 29 ± 1 ) °C  
Relative Humidity ( 66 ± 5 ) %  
Line Voltage ( 220 ± 5 ) Volt

### Condition of this results of Calibration:

- This instrument was calibrated by insert 9 standard thermometer into its chamber and calibration according to W-TE-014 Based on TLAS G-20-1/02-08 (E): Guidelines for Calibration and Checks of Temperature Controlled Enclosures.  
- The temperature scale used was based on ITS - 90.  
- All data show below were final values and the initial data may be obtained upon request.
- Reference Standard Instrument :

Instrument	Model	Serial No./ID No.	Certificate No.	Due Date	Through
Digital Thermometer with sensor	34972A	MY57003188	TE 650469-01	11 June 2023	NATIONAL FOOD INSTITUTE
	RTD	CH#201-209/ RTD#201-209			

- This certificate is traceable to International System of Units (SI Units).
- This certificate was certified only for the instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.
- Condition of Calibrated item : Good

### UUC Description :

Time of Record 1 Hour 9 Minute At 20.0 °C  
Fresh air Damper - Open Position -  
- X Close Fan -  
- Not Available

- Result of Calibration : ☒ Without adjustment ☐ After adjustment

  
3 Aug. 2022

F-CS-012 Revision: 01 Date: 20-04-65



## Calibration Report

**Certificate No.:** 2203876-002-01  
**Equipment:** CHAMBER (Incubator)  
Model: ICP 400 Serial No.: K406.0004  
Resolution: 0.1 °C ID No.: N/A  
Manufacturer: MEMMERT


**Date of Calibration:** 1 August 2022

Page 3 of 3

**Calibration point:** 20.0 °C

**Calibration result:**

Calibration Condition	Temperature (°C)	Relative Humidity (%)	Line Voltage (Volt)
MIN	27.6	61	215.0
MAX	29.5	71	225.0

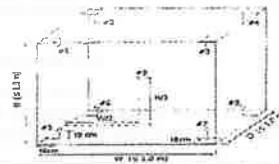

**Table 1 : Reporting of Temperature**

Calibration point (°C)	Measured Temperature (°C) @ Sensor No. (Sensor No.9 is REF)									Uncertainty ± (°C)
	# 1	# 2	# 3	# 4	# 5	# 6	# 7	# 8	# 9	
20.0	20.01	20.09	20.11	20.07	20.18	20.09	20.05	19.99	20.09	0.27

**Table 2 : Reporting of Characterization Result**

UUC* Setting (°C)	UUC* reading (°C)			Stability ± (°C)	Uniformity (°C)	Overall Variation (°C)
	MIN	MAX	Average			
20.0	20.0	20.0	20.0	0.10	0.10	0.37

**Note** The quoted uncertainty include " Stability " and " Loading effect (20% of Temp Uniformity) "

UUC\* = Unit Under Calibration

Stability = One-half of the greatest maximum difference of measured temperatures at any one sensors, for at least half an hour after reaching steady state.

Uniformity = The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time.

Overall Variation = The difference of the maximum and minimum measured temperatures throughout observation time.

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k= 2, providing a level of confidence of approximately 95 %.

----- End -----

FC5-012 Revision: 01 Date: 20-04-65



## Agilent CrossLab Start Up Services Agilent 7890 Gas Chromatograph Preventive Maintenance Checklist

Agilent Preventive Maintenance provides factory recommended service for your analytical instruments to assure reliable operation and the accuracy of your results.

Delivered by highly trained and certified service engineers using genuine Agilent parts and supplies, Agilent Preventive Maintenance provides everything you need to reduce unplanned downtime and keep your systems operating at their peak. This checklist will be completed at the end of the service and provided to you as a record of the preventive maintenance activities.

## Introduction

### Customer Information

- Customers should provide all necessary operating supplies upon request of the engineer.
- A customer representative should be available to the engineer while performing the preventive maintenance procedures.
- Any parts, not included in the Parts Lists section of this document, are not part of the recommended Preventive Maintenance service, nor are they included in the price of this service.
- If a system requires the use of extra or special procedures and/or parts for the maintenance service, then these must be ordered separately and charged as a repair, which may incur additional costs.

### Important Customer Web Links

- For more information about **Agilent Technologies services**, please visit our website using the following URL: <http://www.agilent.com/en-us/products/crosslab-instrument-services/service-repair>
- The **Agilent Community** is an excellent place to get answers, collaborate with others about applications and Agilent products, and find in-depth documents and videos relevant to Agilent technologies. Visit <https://community.agilent.com/welcome>.
- To access **Agilent University**, visit <http://www.agilent.com/crosslab/university/> to learn about training options, which include online, classroom and onsite delivery. A training specialist can work directly with you to help determine your best options.
- A useful **Agilent Resource Center** web page is available, which includes short videos on maintenance, quick lists of consumables for new instruments, and other valuable information. Check out the Resource Page here: <https://www.agilent.com/en-us/agilentresources>.
- Need technical support, FAQs, supplies? – visit our **Support Home page** <http://www.agilent.com/search/support>.
- **Videos** about specific preparation requirements for your instrument can be found by searching the **Agilent YouTube channel** at <https://www.youtube.com/user/agilent>.
- **7890B Manuals** are also available on Agilent.com:
  - **Safety**  
[https://www.agilent.com/cs/library/usermanuals/public/7890B\\_Safety.pdf](https://www.agilent.com/cs/library/usermanuals/public/7890B_Safety.pdf)
  - **Installation and First Startup**  
[https://www.agilent.com/cs/library/usermanuals/Public/7890B\\_Installation.pdf](https://www.agilent.com/cs/library/usermanuals/Public/7890B_Installation.pdf)
  - **Operation Manual**  
[https://www.agilent.com/cs/library/usermanuals/Public/7890B\\_Operation.pdf](https://www.agilent.com/cs/library/usermanuals/Public/7890B_Operation.pdf)
  - **Maintaining Your GC**  
[https://www.agilent.com/cs/library/usermanuals/public/G3430-90052%207890B\\_Maintaining%20Guide.pdf](https://www.agilent.com/cs/library/usermanuals/public/G3430-90052%207890B_Maintaining%20Guide.pdf)

## Service Engineer's Responsibilities

- Contact the customer and ensure that all necessary supplies are available before the preventive maintenance visit.
- Only select those pages that relate to the system or module being serviced.
- Complete empty fields with the relevant information.
- Complete the relevant checkboxes in the checklist using either a "X" or tick mark "✓".
- Check "**Section not applicable**" check boxes to indicate services/tasks not delivered, as appropriate.
- Complete the Preventive Maintenance service in the order of the tasks listed.
- Complete the Service Review section together with the customer.
- Complete the fields for page numbers at the foot of each selected page
- Complete the total number of pages field in the Service Completion section
- **Ask the customer to sign the Service Completion section including the customer's and your signature.**

### Additional Instruction Notes

- Check for any active service notes for this unit. If there are any applicable "Safety" or "Modification Recommended" Service notes, plan to implement the changes on this unit before doing any qualification service.
- Do not implement firmware updates, unless you get approval from the customer and are sure that they are compatible with the instrument control software.

## System Information

- ☒ Check this box if an instrument configuration report is attached instead of completing the table below.

Instrument System Name and ID	6C72908
Instrument System Site and Location	SECOT / INST. 2

List System Component Product Numbers	List the Serial Numbers of each Component
1. 634448	CN 15343147
2. 64513A	CN 11350133
3. 64514A	CN 13020006
4.	
5.	
6.	
7.	
8.	
9.	
10.	

## Preparation

- ☒ Discuss any specific issues with the customer before starting.
- ☒ Review the instrument logbook for recorded problems and comments.
- ☒ Save instrument control settings before starting the procedure.
- ☒ Perform a general inspection of the system for cleanliness.
- ☒ Check for proper installation of parts, assemblies, sensors etc.
- ☒ Check system for required installation of components, settings as defined by current Service Notes.
- ☒ Check for required firmware updates and verify with customers if they would like them installed.
- ☒ Before starting the following procedures, record the Detector Signal Output(s) in the results table. If the GC is turned OFF or in a service mode, comparing the detector outputs before and after the service is not possible.

## Preventive Maintenance Procedure

### Clean and inspect GC

- ☒ Unplug power cord from the power source.
- ☒ Open GC covers and vacuum/remove any dust/debris. Pay particular attention to cooling fans.
- ☒ Inspect internal connectors for proper contact and placement.
- ☒ Reconnect Power to the GC. Power the GC on and verify the power on self-test passed.
- ☒ Verify oven motor spins freely and turns on with the oven door closed; off when the door is opened.
- ☒ Verify operation of all other fans - the inlet and EPC cooling fans.
- ☒ Verify oven intake/outlet flap assembly is operating smoothly while heating and cooling the oven

### Inlet and detector consumable replacement

- ☒ For the inlets installed, perform inlet maintenance as defined in the 7890 manual – "Maintaining Your GC" - for the inlet(s) installed.
- ☒ Replace the split vent trap cartridge filter on units with these inlets: Split/Splitless Capillary (SSL), Multi-Mode Inlet (MMI), Programmed Temperature Vaporizer (PTV), Volatiles Interface (VI).
- ☒ If the inlet system is used in Split Mode with viscous samples, inspect and clean the split vent tube on the inlet and flush or replace the tubing between the inlet and the split vent trap.
- ☒ If the GC includes a Flame Ionization Detector (FID), replace the jet. If the ignitor shows any buildup of sample or corrosion, replace the ignitor. Examine the FID collector and castle assemblies for contamination – clean as necessary.

### Zero Sensors and Leak test

- ☒ Zero all pressure sensors per the procedure in the 7890 "Advanced User Guide".
- ☒ Perform inlet pressure decay test(s) as defined in the 7890 "Troubleshooting Manual". If the PM is done in preparation for an Operational Qualification, then the pressure decay test defined within that protocol can be used for the PM.
- ☒ Record if test passed or failed in the results table.

## ALS Maintenance

- ☐ Section NOT applicable
- ☒ Check all cabling and configuration settings between GC, tray, and injectors.
- ☒ Vacuum or remove any dust, especially around fans.
- ☒ Check operation of all fans.
- ☒ Check syringe for smooth plunger operation.
- ☒ Check for smooth operation of the needle support – clean if necessary

## Restore Instrument

- ☒ Restore the normal operating conditions or customer method using the Data System.
- ☒ Purge the system with carrier flow for 15 minutes
- ☒ Bake out the system, then restore the normal operating conditions
- ☒ After equilibration, check and record the post PM detector signal output values.  
Results should be similar or lower than the detector outputs recorded prior to PM.
- ☒ Perform a chemical checkout. If this is a routine PM, inject the customer's sample using the ALS if applicable. This will act as a final checkout of both the ALS and the GC.

Note: If the PM Service is performed prior to a qualification service, then use the qualification procedure as a guide for final instrument set up and checkout.

## Signature Page

### Service Review

- ☒ Attach available reports/printouts of all tests to this documentation.
- ☒ Record the Preventive Maintenance service activity in the customer's records/logbook.
- ☒ Update/reset instrument maintenance counters as appropriate.
- ☒ Affix the PM sticker to the system or instrument logbook based on the customer's request.
- ☒ Complete the Service Engineer Comments section if there are additional comments.
- ☒ Review with the customer this service, parts replaced, and test results obtained.
- ☐ If the instrument firmware was updated, record the details of the change in the Service Engineer's Comments box or if necessary, in the customer's IQ records.
- ☐ Supply the customer with a copy of the Smart Alerts flyer.
- ☐ Describe Smart Alerts to the customer.
- ☐ Install Smart Alerts if requested.

### 7890 GC Test Results Table

Detector Signal Outputs	Before PM Service	After PM Service
Front detector output	N/A	275.4
Back detector output	N/A	17.6
AUX detector output	N/A	N/A
Pressure decay test	Expected test result	Actual test result
Front inlet pressure decay test	Pass	Pass
Back inlet pressure decay test	Pass	Pass

## 7890 Parts List Table

The following kits are recommended for capillary and purged packed inlets. If this is a general PM and the customer has a preferred set of consumables, you may use the customer's consumables.

Part description	Part number	Product or model# where used	Quantity consumed
SSL Capillary Inlet PM kit, Splitless	5188-6497	7890A/B	-
SSL Capillary Inlet PM kit, split	5188-6496	7890A/B	2
SSL Capillary Ultra Inert Inlet Gold Seal with Washer	5190-6144	7890A/B	-
SSL Capillary Ultra Inert Inlet Splitless Liner - Single taper with Glass Wool	5190-2293	7890A/B	-
SSL Capillary Ultra Inert Inlet Low Pressure Drop Split Liner - with Glass Wool	5190-2295	7890A/B	-
PP Inlet PM kit	5188-6498	7890A/B	-
Split vent trap PM kit, single cartridge (for MMI, PTV & VI)	5188-6495	7890A/B	-
MMI Cleaning Kit	G3510-60820	7890A/B	-
PTV Septumless Head Rebuild Kit	5182-9747	7890A/B	-
PTV Septumless Head Teflon Guide	5182-9748	7890A/B	-
Ignitor (glow plug) assembly with O-ring	19231-60680	7890A/B	1
FID Collector Rebuild/Cleaning Kit	G1531-67000	7890A/B	-
Standard .011-inch FID Jet for capillary FID base	G1531-80560	7890A/B	1
High Temperature .018-inch FID Jet for capillary FID base	G1531-80620	7890A/B	-
Standard .018-inch FID Jet for packed column with packed FID base	18710-20119	7890A/B	-
Standard .011-inch FID Jet for capillary column with packed/adaptable FID base	19244-80560	7890A/B	-
High Temperature .018-inch FID Jet for capillary column with packed/adaptable FID base	19244-80620	7890A/B	-
NPD Jet, universal fit, .011-inch ID	G1534-80580	7890A/B	-
NPD Jet, universal fit, .011-inch ID Extended tip	G1534-80590	7890A/B	-
SSL Capillary Ultra Inert Inlet Gold Seal with Washer	5190-6144	7890A/B	-
SSL Capillary Ultra Inert Inlet Splitless Liner - Single taper with Glass Wool	5190-2293	7890A/B	-
**FID Collector Replacement Kit, if needed	G1531-67001	7890A/B	-

## Service Engineer Comments

If there are any specific points you wish to note as part of performing the service or other items of interest for the customer, please write include them in this box.

NA

## Service Completion

Service request number 6005212573

Date service completed 12-June-2022

Agilent signature Sam Pop T.

Customer signature Mamm. Samanipale

Total number of pages in this document 9



## Agilent CrossLab Start Up Services

# Agilent GCMS Preventive Maintenance Checklist



Agilent Preventive Maintenance provides factory recommended service for your analytical instruments to assure reliable operation and the accuracy of your results.

Delivered by highly trained and certified service engineers using genuine Agilent parts and supplies, Agilent Preventive Maintenance provides everything you need to reduce unplanned downtime and keep your systems operating at their peak. This checklist will be completed at the end of the service and provided to you as a record of the preventive maintenance activities.

## Introduction

Select the appropriate PM to be done and then perform the checklist under that section

- ☐ Interim Preventive Maintenance 6 months
- ☒ Major Preventive Maintenance Yearly

This checklist covers the following model(s):

Type	Model
SQ	5973 Series MSD
SQ	5975 Series MSD
SQ	5977 Series MSD
TQ	7000 Series MS/MS
TQ	7010 Series MS/MS
QTOF	7200 Series QTOF
QTOF	7250 Series QTOF

## Customer Information

- Customers should provide all necessary operating supplies upon request of the engineer.
- A customer representative should be available to the engineer while performing the preventive maintenance procedures. Customers are responsible for regular maintenance and are encouraged to observe the service representative.
- Any parts not included in the Parts Lists section of this document are not part of the recommended Preventive Maintenance service nor are they included in the price of this service.
- If a system requires the use of extra or special procedures and/or parts for the maintenance service, then these must be ordered separately and charged as a repair, which may incur additional costs.

### Important Customer Web Links

- For more information about Agilent Technologies services, please visit our website using the following URL: <http://www.agilent.com/en-us/products/crosslab-instrument-services/service-repair>
- To access Agilent University, visit <http://www.agilent.com/crosslab/university/> to learn about training options, which include online, classroom and onsite delivery. A training specialist can work directly with you to help determine your best options.
- A useful Agilent Resource Center web page is available, which includes short videos on maintenance, quick lists of consumables for new instruments, and other valuable information. Check out the Resource Page here: <https://www.agilent.com/en-us/agilentresources>
- Need technical support, FAQs, supplies? – visit our Support Home page at <http://www.agilent.com/search/support>
- Get answers. Share insights. Build connections: Join the Agilent Community at <https://community.agilent.com/welcome>

### Service Engineer's Responsibilities

- Contact the customer and ensure that all necessary supplies are available before the preventive maintenance visit.
- Complete empty fields with the relevant information.
- Complete the relevant checkboxes in the checklist using either a "X" or tick mark "✓".
- Check "Section not applicable" check boxes to indicate services/tasks not delivered, as appropriate.
- Complete the Preventive Maintenance services in the most logical order relevant to the individual system service in the order of the tasks listed.
- Complete the Service Review section together with the customer.
- Ask the customer to sign the Service Completion section including the customer's and your signature.

### Additional Instruction Notes

- Preventive maintenance is a factory recommended procedure designed to reduce the likelihood of electromechanical failures. Failure to perform preventive maintenance may reduce the long-term reliability of certain instruments and systems. Two preventative maintenances (PMs) per year are recommended, the Major PM Service will be performed annually with an Interim PM performed 6 months after the Major PM.

### System Information

- ☐ Check this box if an instrument configuration report is attached instead of completing the table.

Instrument System Name and ID 5975C MSD / GCMS-FID  
Instrument System Site and Location SECOT Co., Ltd.

### List System Component Product Numbers List the Serial Numbers of each Component

1. G3172A	US74838080
2.	
3.	
4.	
5.	
6.	
7.	
8.	

### Preparation

- ☒ Discuss any specific issues with the customer before starting.
- ☒ Review the instrument logbook for recorded problems and comments.
- ☒ Save instrument control settings before starting the procedure.
- ☒ Perform a general inspection of the system for cleanliness.
- ☒ Check for proper installation of parts, assemblies, sensors etc.
- ☒ Check system for required installation of components and settings as defined by current Service Notes
- ☒ Check for firmware updates and verify with customers if they would like them installed. Firmware update(s) are strongly recommended.

### Customer Responsibilities

Customers should ensure that all necessary operating supplies, consumables, and usage-dependent items such as gases, vials, syringes, calibrant solution and solvents required for successful preventive maintenance are available. A customer representative should be available while the preventive maintenance is being performed.

## Important notice for customers

The customer should complete the following before the Support Provider arrives on site:

- ☒ Perform an autotune and retain the printed tune report just prior to the start of the PM to verify performance of the equipment.

Note: it is recommended to have the customer run the autotune and tune evaluation prior to the PM and then start the vent cycle so that the instrument will be ready for the service representative.

## Definition of the Task/Recommended items within the document

Task	Recommended	
Yes	No	Interim / Major / As needed
<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>

Yes selected means that the task was done or the part was required.

No selected means that the task was not done or the part was not required.

Interim selected means that this task is recommended to be done at 6-month intervals.

Major selected means that this task is recommended to be done yearly; if the customer would like a service to be done at the 6-month interval then the service could be purchased.

As needed selected means that the task was done or the part was used as needed. For example, there could be two types of filters that could be used and this was the one selected.

## Preventive Maintenance Procedures

Yes/No	Interim/Major	Description
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Perform general inspection of system for cleanliness
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Discuss any problems the customer is having with the instrument
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Review customer maintenance records and exclude maintenance on recently serviced items
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Review the most recent autotune report. This will give a starting point for evaluating spectral peaks, baseline noise, peak shape, mass assignments and resolution

Yes/No	Interim/Major	GCMS
Yes/No	Interim/Major	Description
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Record instrument model no.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Record instrument serial no.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Record Rough Vacuum
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Record Manifold Vacuum
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Type of Column installed

Yes/No	Interim/Major	System Checks
Yes/No	Interim/Major	Description
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Verify that calibration peaks were seen prior to starting the PM
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Vent the instrument
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Inspect vacuum hoses, pump, exhaust tubing, and power cords for excessive wear.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Visually inspect calibrant levels – PFTBA PFDTD (if appl.), IRM (if appl.). Refill if available.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Look for any obvious external damage or problems.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Clean air intake(s). Cosmetic cover(s) may need to be removed
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Verify system line voltage meets instrument specifications: Yes <input checked="" type="checkbox"/> No <input type="checkbox"/>

Yes/No	Interim/Major	Wet Mechanical vacuum pumps
Yes/No	Interim/Major	Description
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Check for evidence of oil leakage. Check pump gasket for leakage.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Drain and replace mechanical pump oil.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Replace Oil Mist Filter if applicable.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Discuss with customer the need for more frequent oil changes if the oil is dirty
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Don't use mist filters with Chemical Ionization
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Perform anti-suckback valve test. Power on until side plate is held closed, power off and check that side plate holds closed. Visually confirm that no oil returns up vacuum hose.
Yes/No	Interim/Major	Dry Mechanical vacuum pumps - Diaphragm
Yes/No	Interim/Major	Description
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Check for evidence of poor vacuum – Turbo power demand, poor manifold vacuum, etc.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Clear air flow paths of dust.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	If vacuum is poor, then replace the diaphragm pump
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Perform anti-suckback valve test. Power on until side plate is held closed, power off and check that side plate holds closed.

Yes/No	Interim/Major	Dry Mechanical vacuum pumps - Scroll
Yes/No	Interim/Major	Description
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Replace the tips seal on the IDP pump.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Check for evidence of poor vacuum – Turbo power demand, poor manifold vacuum, etc.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Replace the Exhaust Filter if required.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Discuss with customer the need for more frequent changes, if needed.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Inform customer that pump gas ballast should be installed all the time.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Perform anti-suckback valve test. Power on until side plate is held closed, power off and check that side plate holds closed.

Yes/No	Interim/Major	Cleaning System and Filters
Yes/No	Interim/Major	Description
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Fans
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Remove dust from fans and vent covers.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Verify fans are functional and that there is enough space around the instrument for proper cooling.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Source cleaning
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Open analyzer and remove the source.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Disassemble, Clean, Re-assemble source.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Re-install source and close analyzer.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Filters
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Replace RMSH-2 Helium gas filter – if applicable.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Replace RMSN-2 Nitrogen gas filter – if applicable.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	Replace RMSHY-2 Hydrogen gas filter – if applicable.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	CP17988 – Gas Clean Carrier Gas Kit for 7890 for Nitrogen or Helium, Bracket, Mount, and Filter – if applicable.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	CP17974 – Gas Clean Filter Kit GC/MS 1/8"; Mount and Filter – if applicable.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	CP17973 – Gas Clean Filter; Replacement Filter – if applicable.
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	S190-9071 – Methane Gas Filter – if applicable.

Guidance: if gas filter is replaced, write the change date on the filter using a permanent marker.

**Guidance:** If the PM Service is performed prior to a qualification service, then use the qualification procedure as a guide for final instrument setup and checkout.

- ☒ Attach available reports/printouts of all tests to this documentation.
- ☒ Record the Preventive Maintenance service activity in the customer's records/logbook. Record the PM event in the Smart Alerts logbook, if applicable.
- ☒ Update/reset instrument maintenance counters as appropriate.
- ☒ Affix the PM sticker to the system or instrument logbook based on the customer's request.
- ☒ Complete the Service Engineer Comments section if there are additional comments.
- ☒ Review this service, parts replaced, and test results obtained with the customer.
- ☐ If the instrument firmware was updated, record the details of the change in the Service Engineer's Comment box. Systems in a compliant environment may need additional documentation.

[illegible]☐ Section not applicable[illegible]

## Signature Page

### Service Engineer Comments (optional)

N/A

### Service Completion

Service request number 6005218299 Date service completed 17 JUN, 2022

Agilent signature Shirya T. Customer signature Shirwan C.

Total number of pages in this document 12

## Parts – As needed as part of the PM

### Common MS Filters and Seals – 5973/5975/5977/7000/7010/7200/7250 Series

Yes/No	Interim/Major/As needed	Supplies	Part number
<input type="checkbox"/>	<input type="checkbox"/>	Helium gas filter – if required	RMSH-2
<input type="checkbox"/>	<input type="checkbox"/>	Nitrogen gas filter – if required	RMSN-2
<input type="checkbox"/>	<input type="checkbox"/>	Big Universal Trap, 1/8" fittings, Hydrogen, if required	RMSHY-2
<input type="checkbox"/>	<input type="checkbox"/>	Gas Clean Carrier Gas Kit for 7890 for Nitrogen or Helium; Bracket, Mount and Filter – if required	CP17988
<input type="checkbox"/>	<input type="checkbox"/>	Gas Clean Filter Kit GC/MS 1/8 in (complete replacement kit) – if required	CP17974
<input type="checkbox"/>	<input type="checkbox"/>	Gas Clean GS/MS Filter – if required	CP17973
<input type="checkbox"/>	<input type="checkbox"/>	Chemical Ionization Gas Purifier (CI systems) – if required	5190-9071
<input checked="" type="checkbox"/>	<input type="checkbox"/>	Agilent AVF Platinum, 1 quart	5191-5851

Gas filters need to be changed only if required

### MS Maintenance Supplies for 5973/5975/5977 Series

Yes/No	Interim/Major/As needed	Supplies	Part number
<input type="checkbox"/>	<input type="checkbox"/>	Diffusion pump fluid (Diffusion Pump Models)	6040-0809 Qty 2
<input type="checkbox"/>	<input type="checkbox"/>	IDP-3 Tip Seal Replacement Kit (IDP-3 Dry Pump Models)	G7077-67018
<input type="checkbox"/>	<input type="checkbox"/>	IDP-3 Tip Seal Replacement Kit (no tools – CSD P/N)	5190-9561
<input type="checkbox"/>	<input type="checkbox"/>	IDP-3 Tip Seal Replacement Kit (no tools – VPD P/N)	IDP3TS
<input type="checkbox"/>	<input type="checkbox"/>	Filter element for IDP-3	REPLSLRFILTER2
<input type="checkbox"/>	<input type="checkbox"/>	DS42 Oil Mist Eliminator 3/4G & 3/8	SR03706556
<input type="checkbox"/>	<input type="checkbox"/>	Exhaust oil mist trap (thread) Edwards/Pfeiffer	G1099-80039

### MS Maintenance Supplies for 7000/7010 Series

Yes/No	Interim/Major/As needed	Supplies	Part number
<input type="checkbox"/>	<input type="checkbox"/>	Nitrogen gas filter	RMSN-2
<input type="checkbox"/>	<input type="checkbox"/>	IDP-10 Tip Seal Replacement Kit (IDP-10 Dry Scroll Pump Models)	G7004-67023
<input type="checkbox"/>	<input type="checkbox"/>	IDP-10 Tip Seal Replacement Kit (no tools – VPD P/N)	X3807-67000
<input type="checkbox"/>	<input type="checkbox"/>	Oil Mist Filter RV5	G6600-80043
<input type="checkbox"/>	<input type="checkbox"/>	Filter element for the IDP-10	REPLSLRFILTER1

### MS Maintenance Supplies for 7200/7250 Series

Yes/No	Interim/Major/As needed	Supplies	Part number
<input type="checkbox"/>	<input type="checkbox"/>	Nitrogen gas filter – if required	RMSN-2
<input type="checkbox"/>	<input type="checkbox"/>	RIS Probe Maintenance Kit (7200 Series only)	G7005-60170
<input type="checkbox"/>	<input type="checkbox"/>	DS202 Oil Mist Eliminator	SR03706800
<input type="checkbox"/>	<input type="checkbox"/>	IDP-15 Tip Seal Replacement Kit (IDP-15 Dry Pump Models)	5190-9613
<input type="checkbox"/>	<input type="checkbox"/>	IDP-15 Tip Seal Replacement Kit (no tools – VPD P/N)	X3815-67000
<input type="checkbox"/>	<input type="checkbox"/>	Filter element, for SH-110/SH-112/IDP-15 exhaust silencer	REPLSLRFILTER
<input type="checkbox"/>	<input type="checkbox"/>	DS 3/8 MAG. PLUG AND GASKET	SR03701824

### MS Maintenance Supplies for JetClean

Yes/No	Interim/Major/As needed	Supplies	Part number
<input type="checkbox"/>	<input type="checkbox"/>		

<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>	Big Universal Trap, 1/8" fittings, Hydrogen, if required	RMSHY-2
--------------------------	--------------------------	--------------------------	-------------------------------------	--	---------

## Consumable Parts Reference – Purchasable by customer, not included as part of PM

### Common MSD Maintenance Supplies 5973/5975/5977/7000/7010/7200/7250 Series

Yes/No	Interim/Major/As needed	Description	Part number
<input checked="" type="checkbox"/>	<input type="checkbox"/>	El High Temperature Filaments	G7005-60061 Qty 2
<input type="checkbox"/>	<input type="checkbox"/>	HES El Filaments	G7002-60001
<input type="checkbox"/>	<input type="checkbox"/>	LE-El Filaments	G3850-60021
<input type="checkbox"/>	<input type="checkbox"/>	Cl High Temperature Filament – all MSDs	G7005-60072
<input type="checkbox"/>	<input type="checkbox"/>	PFTBA GCMS Tuning Standard calibrant	05971-60571
<input type="checkbox"/>	<input type="checkbox"/>	PFTD calibrant, 1 mL	8500-8510
<input type="checkbox"/>	<input type="checkbox"/>	PFET, IRM calibrant for GC QTOF 0.5 mL	5190-0531

### MSD Maintenance Supplies 5973/5975/5977 Series

Yes/No	Interim/Major/As needed	Description	Part number
<input type="checkbox"/>	<input type="checkbox"/>	Cl Interface tip seal (tip and spring combo)	G1999-60412
<input type="checkbox"/>	<input type="checkbox"/>	Cl Interface tip seal (tip only)	G3870-20542
<input type="checkbox"/>	<input type="checkbox"/>	Cl Interface tip seal spring (spring only)	G1999-20023
<input type="checkbox"/>	<input type="checkbox"/>	Repeller insulator	G1099-20133 Qty 2
<input type="checkbox"/>	<input type="checkbox"/>	Lens insulator/holder (HES)	G7002-20074
<input type="checkbox"/>	<input type="checkbox"/>	Ring heater/sensor assembly (HES)	G7002-60043
<input type="checkbox"/>	<input type="checkbox"/>	Ceramic insulator for Extractor (HES)	G7002-20064
<input type="checkbox"/>	<input type="checkbox"/>	Transfer-Line Tip Cap, Threaded	G3870-20547
<input type="checkbox"/>	<input type="checkbox"/>	Transfer-Line Tip Base, Threaded	G3870-20548

### MS Maintenance Supplies for 7000/7010 Series

Yes/No	Interim/Major/As needed	Description	Part number
<input type="checkbox"/>	<input type="checkbox"/>	Cl Interface tip seal - 7000	G1999-60412
<input type="checkbox"/>	<input type="checkbox"/>	Cl Interface tip seal - 7010	G7002-60412
<input type="checkbox"/>	<input type="checkbox"/>	Cl Interface tip seal (tip only)	G3870-20542
<input type="checkbox"/>	<input type="checkbox"/>	Cl Interface tip seal spring (spring only)	G1999-20023
<input type="checkbox"/>	<input type="checkbox"/>	Repeller insulator - 7000	G1099-20133 Qty 2
<input type="checkbox"/>	<input type="checkbox"/>	Lens insulator/holder (HES)	G7002-20074
<input type="checkbox"/>	<input type="checkbox"/>	Ring heater/sensor assembly (HES)	G7002-60043
<input type="checkbox"/>	<input type="checkbox"/>	Ceramic insulator for Extractor (HES)	G7002-20064
<input type="checkbox"/>	<input type="checkbox"/>	Transfer-Line Tip Cap, Threaded	G3870-20547
<input type="checkbox"/>	<input type="checkbox"/>	Transfer-Line Tip Base, Threaded	G3870-20548

### MS Maintenance Supplies for 7200 Series

Yes/No	Interim/Major/As needed	Description	Part number
<input type="checkbox"/>	<input type="checkbox"/>	Extractor Lens Insulator	G7005-20133
<input type="checkbox"/>	<input type="checkbox"/>	Ion Focus Insulator	G7005-20442
<input type="checkbox"/>	<input type="checkbox"/>	Ring Heater/Sensor Assembly	G7005-60110
<input type="checkbox"/>	<input type="checkbox"/>	RIS Xfer Tip	G7005-20542
<input type="checkbox"/>	<input type="checkbox"/>	RIS Xfer Tip Spring	G7005-20024

### MS Maintenance Supplies for 7250 Series

Yes/No	Interim/Major/As needed	Description	Part number
<input type="checkbox"/>	<input type="checkbox"/>	Lens insulator/holder (HES)	G7002-20074
<input type="checkbox"/>	<input type="checkbox"/>	Ring heater/sensor assembly (HES)	G7002-60043
<input type="checkbox"/>	<input type="checkbox"/>	Ceramic insulator for Extractor (HES)	G7002-20064
<input type="checkbox"/>	<input type="checkbox"/>	Transfer-Line Tip Cap, Threaded	G3870-20547
<input type="checkbox"/>	<input type="checkbox"/>	Transfer-Line Tip Base, Threaded	G3870-20548
<input type="checkbox"/>	<input type="checkbox"/>	El Extractor Transfer Tip	G3870-20542
<input type="checkbox"/>	<input type="checkbox"/>	Cl Tip Compression Spring	G1999-20023

### MS Maintenance Supplies for Intuvo 9000 MS Systems

Yes/No	Interim/Major/As needed	Description	Part number
<input type="checkbox"/>	<input type="checkbox"/>	Swaged MS Tail - Packaged	G4590-60009
<input type="checkbox"/>	<input type="checkbox"/>	Swaged MS Tail (HES) - Packaged	G4590-60109

### Common MS Maintenance Supplies

Yes/No	Interim/Major/As needed	Description	Part number
<input type="checkbox"/>	<input type="checkbox"/>	Abrasive paper, 30 um	5061-5896
<input type="checkbox"/>	<input type="checkbox"/>	Alumina powder	393706201
<input type="checkbox"/>	<input type="checkbox"/>	Cloths, clean (pkg of 15)	05980-60051
<input type="checkbox"/>	<input type="checkbox"/>	Cloths, cleaning (pkg of 300)	9310-4828
<input type="checkbox"/>	<input type="checkbox"/>	Cotton swabs (pkg of 100)	5080-5400
<input type="checkbox"/>	<input type="checkbox"/>	Gloves, clean, large	8650-0030
<input type="checkbox"/>	<input type="checkbox"/>	Gloves, clean, small	8650-0029