

ภาคผนวก จ

ใบแสดงการตรวจเทียบเครื่องมือ

CERTIFICATE OF ANALYSIS

Grade of Product: EPA Protocol

Part Number: E04NI99E15A0KDC	Reference Number: 82-124380327-1
Cylinder Number: CC414868	Cylinder Volume: 144.4 Cubic Feet
Laboratory: ASG - Riverton - NJ	Cylinder Pressure: 2015 PSIG
PGVP Number: B52013	Valve Outlet: 660
Gas Code: CO,NO,SO2,BALN	Certification Date: Jul 15, 2013

Expiration Date: Jul 15, 2021

Certification performed in accordance with "EPA Traceability Protocol for Assay and Certification of Gaseous Calibration Standards (May 2012)" document EPA 600/R-12/531, using the assay procedures listed. Analytical Methodology does not require correction for analytical interference. This cylinder has a total analytical uncertainty as stated below with a confidence level of 95%. There are no significant impurities which affect the use of this calibration mixture. All concentrations are on a volume/volume basis unless otherwise noted.

Do Not Use This Cylinder below 100 psig, i.e. 0.7 megapascals.

ANALYTICAL RESULTS

Component	Requested Concentration	Actual Concentration	Protocol Method	Total Relative Uncertainty	Assay Dates
NOX	50.00 PPM	52.49 PPM	G1	+/- 1% NIST Traceable	07/05/2013, 07/15/2013
NITRIC OXIDE	50.00 PPM	52.46 PPM	G1	+/- 1% NIST Traceable	07/05/2013, 07/15/2013
SULFUR DIOXIDE	50.00 PPM	50.64 PPM	G1	+/- 1.0% NIST Traceable	07/05/2013, 07/15/2013
CARBON MONOXIDE	0.5000 %	0.4984 %	G1	+/- 0.4% NIST Traceable	07/07/2013
NITROGEN	Balance				

CALIBRATION STANDARDS

Type	Lot ID	Cylinder No	Concentration	Uncertainty	Expiration Date
NTRM	13060232	CC401984	4950 PPM CARBON MONOXIDE/NITROGEN	+/- 0.4%	Feb 15, 2019
PRM	12312	680179	10.01 PPM NITROGEN DIOXIDE/NITROGEN	+/- 2.0%	Feb 14, 2012
NTRM	12060813	CC281093	49.95 PPM NITRIC OXIDE/NITROGEN	+/- 0.8%	Dec 16, 2017
GMIS	124206889106	CC322664	4.879 PPM NITROGEN DIOXIDE/NITROGEN	+/- 2.0%	Apr 08, 2016
NTRM	12061832	CC352180	50.10 PPM SULFUR DIOXIDE/NITROGEN	+/- 1.0%	Apr 24, 2018

The SRM or PRM noted above is only in reference to the GMIS used in the assay and not part of the analysis.

ANALYTICAL EQUIPMENT

Instrument/Make/Model	Analytical Principle	Last Multipoint Calibration
Siemens Ultramat 6 N1C8180 COHIGH	NDIR	Jul 08, 2013
Nicolet 6700 APW1100391 NO	FTIR	Jun 24, 2013
Nicolet 6700 APW1100391 NO2	FTIR	Jun 24, 2013
Nicolet 6700 APW1100391 SO2	FTIR	Jul 10, 2013

Triad Data Available Upon Request

Notes:

C. Modyelewski

Approved for Release



NO Analyzer Performance Test

Date : 8 Feb 21

Temp: (°C) 25

Barometric Pressure: Pb (mmHg) 760

Analyzer Type :	No
Brand :	API
Model :	200A
S/N :	1523

Dilutor	: Teledyne 700E 587
Zero Air	: M701 S/N 1044
STD GAS	: EB0108319

NO Single Point Calibration

Supply Gas	Ref Value	Analyzer Disp.	Zero-Span Error %	Slope - Offset
Zero	0.0	0.4	-	-
Span	450.0	445.0	-1.11	0.993

NO MultiPoint Calibration

Ref Value	Analyzer Disp.	Output Difference		
		Diff	Percent Diff	Percent Diff abs.
0.00	0.40	0.4	-	-
100.00	97.50	-2.5	-2.5	2.5
200.00	194.20	-5.8	-2.9	2.9
400.00	397.60	-2.4	-0.6	0.6
			Average Diff (%)	2.0

Calibrated by :

Approved by :



Nox Analyzer Performance Test

Date : 8 Feb 21

Temp: (°C) 25

Barometric Pressure: Pb (mmHg) 760

Analyzer Type :	Nox
Brand :	API
Model :	200A
S/N :	1523

Dilutor	: Teledyne 700E 587
Zero Air	: M701 S/N 1044
STD GAS	: EB0108319

NOx Single Point Calibration

Supply Gas	Ref Value	Analyzer Disp.	Zero-Span Error %	Slope - Offset
Zero	0.0	1.0	-	-
Span	450.0	448.0	-0.44	0.995

NOx MultiPoint Calibration

Ref Value	Analyzer Disp.	Output Difference		
		Diff	Percent Diff	Percent Diff abs.
0.00	1.00	1.0	-	-
100.00	99.10	-0.9	-0.9	0.9
200.00	197.30	-2.7	-1.3	1.3
400.00	399.10	-0.9	-0.2	0.2
			Average Diff (%)	0.8

Calibrated by :

Approved by :



NO Analyzer Performance Test

Date : 11 Jan 21

Temp: (°C) 25

Barometric Pressure: Pb (mmHg) 760

Analyzer Type :	No
Brand :	Teledyne
Model :	T200
S/N :	111

Dilutor	: Teledyne 700E 587
Zero Air	: M701 S/N 1044
STD GAS	: EB0108319

NO Single Point Calibration

Supply Gas	Ref Value	Analyzer Disp.	Zero-Span Error %	Slope - Offset
Zero	0.0	0.1	-	-
Span	450.0	449.7	-0.07	0.999

NO MultiPoint Calibration

Ref Value	Analyzer Disp.	Output Difference		
		Diff	Percent Diff	Percent Diff abs.
0.00	0.10	0.1	-	-
100.00	101.10	1.1	1.1	1.1
200.00	202.00	2.0	1.0	1.0
400.00	402.00	2.0	0.5	0.5
			Average Diff (%)	0.9

Calibrated by :

Approved by :



Nox Analyzer Performance Test

Date : 11 Jan 21

Temp: (°C) 25

Barometric Pressure: Pb (mmHg) 760

Analyzer Type :	Nox
Brand :	Teledyne
Model :	T200
S/N :	111

Dilutor :	Teledyne 700E 587
Zero Air :	M701 S/N 1044
STD GAS :	EB0108319

NOx Single Point Calibration

Supply Gas	Ref Value	Analyzer Disp.	Zero-Span Error %	Slope - Offset
Zero	0.0	0.3	-	-
Span	450.0	453.0	0.67	1.006

NOx MultiPoint Calibration

Ref Value	Analyzer Disp.	Output Difference		
		Diff	Percent Diff	Percent Diff abs.
0.00	0.30	0.3	-	-
100.00	101.30	1.3	1.3	1.3
200.00	203.40	3.4	1.7	1.7
400.00	405.00	5.0	1.3	1.3
			Average Diff (%)	1.4

Calibrated by : Wittaya KApproved by : [Signature] V



CONTROL UNIT CALIBRATION

(Metric units, mm)

Date 8 Jan 21

Initial Final Average
Barometric press, Pb 758 758 758 mmHg

Dry Gas Meter Data

Console No. M50-09

Metering System ID

DGM Number 972135

DGM Model ES-110

Calibrated by : Montri P.

Reference Dry Gas Meter Data

Serial No. 358794

Model S110

Correction factor (Yr) 1.0096

Last Calibration Date 19 Dec 20

Orifice manometer setting, ΔH mm H2O	Ref. DGM Volume V _r Liters	DGM Volume V _m Liters	Temperature (°C)				Time ⊙ min	DGM Correction factor (Y)	ΔH@ mm
			Ref DGM T _r	Dry Gas Meter					
				Inlet T _i	Outlet T _o	Avg T _m			
12.5	100.1	98.6	25	25	25	25.0	8.28	1.0043	40.1199
25.0	100.1	99.4	25	25	25	25.0	6.01	0.9949	42.2999
50.0	100.1	99.7	25	25	25	25.0	4.38	0.9898	44.8705
76.0	100.2	100.8	25	25	26	25.5	3.30	0.9789	38.5321
100.0	100.0	100.3	25	26	26	26.0	3.30	0.9815	46.3993
150.0	100.2	99.4	25	26	26	26.0	2.53	0.9877	44.6649

Average 0.9895 42.8144

Approved by : 
(Miss Katesarin Vorradetwittaya)



PITOT TUBE CALIBRATION

Calibration Location: SECOT

Calibration Date : Dec 19, 20

Calibrated duct No.: 1

Calibration Standard Pitot tube data

Pitot No. : Std-01

Coefficient (Cp) : 1

Type S Pitot No. : PS20-03

Calibrated by : Mr. Montri P.

A Side Calibration

Run No.	ΔP_{std} (mm H ₂ O)	ΔP_s (mm H ₂ O)	Cp(s)	Deviation, δ Cp(s) - Cp(A)
1	7.50	10.50	0.8452	-0.0034
2	7.50	10.25	0.8554	0.0068
3	7.50	10.50	0.8452	-0.0034

C_{P(A),avg} 0.8486

B Side Calibration

Run No.	ΔP_{std} (mm H ₂ O)	ΔP_s (mm H ₂ O)	Cp(s)	Deviation, δ Cp(s) - Cp(B)
1	7.50	10.50	0.8452	0.0000
2	7.50	10.50	0.8452	0.0000
3	7.50	10.50	0.8452	0.0000

C_{P(B),avg} 0.8452

| CP(A)-CP(B) | = 0.0034

C_{P(Avg)} = 0.8469

Approved by :

(Miss Katesarin Vorradetwittaya)

*** δ must be ≤ 0.01 for the test to be acceptable ***
 *** | Cp(A)-Cp(B) | must also be < 0.01 if average of Cp(A) and Cp(B) is to be used ***



CONTROL UNIT CALIBRATION

(Metric units, mm)

Date 8 Jan 21

	Initial	Final	Average	
Barometric press, Pb	758	758	758	mmHg

Dry Gas Meter Data

Console No. M50-08

Metering System ID

DGM Number 971415

DGM Model ES-110

Calibrated by : Montri P.

Reference Dry Gas Meter Data

Serial No. 358794

Model S110

Correction factor (Yr) 1.0096

Last Calibration Date 19 Dec 20

Orifice manometer setting, ΔH mm H2O	Ref. DGM Volume V _r Liters	DGM Volume V _m Liters	Temperature (°C)				Time Θ min	DGM Correction factor (Y)	ΔH@ mm
			Ref DGM T _r	Dry Gas Meter					
				Inlet T _i	Outlet T _o	Avg T _m			
12.5	100.4	100.6	25	25	25	25.0	8.25	1.0058	38.1338
25.0	100.6	99.8	25	25	25	25.0	6.33	1.0154	44.7234
50.0	100.0	99.7	25	25	25	25.0	4.03	1.0075	36.6528
76.0	100.0	98.7	25	26	26	26.0	3.55	1.0191	43.0866
100.0	100.1	100.6	25	25	26	25.5	3.55	0.9964	41.2182
150.0	100.0	99.7	25	25	26	25.5	2.42	0.9995	39.5179

Average	1.0073	40.5554
---------	--------	---------

Approved by :

(Miss Katesarin Vorradetwittaya)



PITOT TUBE CALIBRATION

Calibration Location: SECOT

Calibration Date : Dec 19, 20

Calibrated duct No.: 1

Calibration Standard Pitot tube data

Pitot No. : Std-01

Coefficient (Cp) : 1

Type S Pitot No. : PS20-02

Calibrated by : Mr. Montri P.

A Side Calibration

Run No.	ΔP_{std} (mm H ₂ O)	ΔP_s (mm H ₂ O)	Cp(s)	Deviation, δ Cp(s) - Cp(A)
1	7.50	10.50	0.8452	0.0000
2	7.50	10.50	0.8452	0.0000
3	7.50	10.50	0.8452	0.0000

C_{P(A),avg} 0.8452

B Side Calibration

Run No.	ΔP_{std} (mm H ₂ O)	ΔP_s (mm H ₂ O)	Cp(s)	Deviation, δ Cp(s) - Cp(B)
1	7.50	10.75	0.8353	0.0000
2	7.50	10.75	0.8353	0.0000
3	7.50	10.75	0.8353	0.0000

C_{P(B),avg} 0.8353

| CP(A)-CP(B) | = 0.0099

C_{P(Avg)} = 0.8402

Approved by : 
(Miss Katesarin Vorradetwittaya)

*** δ must be ≤ 0.01 for the test to be acceptable ****** | Cp(A)-Cp(B) | must also be < 0.01 if average of Cp(A) and Cp(B) is not be used ***

THE LINDE GROUP

Linde

Certificate of Analysis

Special Gases Mixture

Customer Details

Name:

LABIX CO.,LTD.

Address:

House number.105/12 Moo 2,Sukhumvit Rd., T.
Thung Suk La A.Sriracha Chon Buri 20230

Customer Tag No.:

Certificate Details

Number:

2237/18

Date of Issue:

28-Jun-2018

Expired date:

27-Jun-2022

Material Details

Production Order:

90149087

Material Code:

478100-J-62

Cylinder No.:

S197614

Gas content:

6.52 M³ (nominal)

Filling pressure:

145 bar (g)

Valve:

CGA 590 BRASS

Cylinder Owner:

LINDE

Cylinder Material:

STEEL

Cylinder Size:

47 L

Laboratory Report

Component	Normal Concentration	Analysis Result ¹	Uncertainty ²	Method of Analysis ³
Oxygen	8.00%	7.74%	± 2% relative	(1) SG-O-01
Nitrogen	Balance			

Recommend usage condition

Minimum utilization: 5% of actual content or before expire date whichever comes first.

Storage condition: Keep in well ventilation and secure area.

Comments

Note:

1. All results expressed in this report are on mole/mole basis, unless otherwise specified.
2. The reported expanded uncertainty is based on a standard uncertainty multiplied by a coverage factor k=2, providing a level of confidence of approximately 95%. The measurement of this material is traceable to the SI through the reference gas standard which is traceable to Swiss National Standard of Mass or other recognised national metrology institutes.
3. (1) Gas Chromatography, (2) Paramagnetic Oxygen Analyzer, (3) Electrochemical Oxygen Analyzer, (4) Electrochemical Moisture Analyzer, (5) Total Hydrocarbon Analyzer, (6) Other - Specified

Sukanya Kamutharat

Signatory for and on behalf of Linde (Thailand) Co., Ltd.

Page 1 of 1

This report shall not be reproduced except in full

PB-002/F004

Iss:H/2, 01 March 2018

บริษัท ลินด์ (ประเทศไทย) จำกัด (มหาชน)

เลขที่ใบอนุญาตประกอบกิจการ 0107537000785

ชั้น 15 อาคารทาวเวอร์ เอ 2/3 หมู่ 14 ถนนบางนา-ตราด กม. 6.5 แขวงบางนา

อ.บางพลี จ.สมุทรปราการ 10540 โทรศัพท์ (66) 2338-6100 โทรสาร (66) 2338-6333

โรงงานเวลโกรว์: 105 หมู่ 5 ต.บางพลีใหญ่ อ.บางพลี จ.สมุทรปราการ 24180

โทรศัพท์ (66) 38.570-479-93

โทรสาร (66) 38.570-323

Linde (Thailand) Public Company Limited

PLC., Registration no.0107537000785

15th Floor, Bangna Tower A, 2/3 Moo 14, Bangna Trad KM. 6.5 Road, Bangkaew

Bangplee, Samutprakarn 10540, Tel (66) 2338-6100 Fax (66) 2338-6333

Wellgrow Plant: 105 Moo 5, T.Bangsamak, A.Bangpakong, Chachoengsao 24180

Thailand, Tel (66) 38.570-479-93

Fax (66) 38.570-323

THE LINDE GROUP

Linde

Certificate Of Analysis

Special Gases Mixture

Customer Details

Name:

Secot Co., Ltd.

Address:

239, Rimklongprapa Rd., Bangsue, Bangkok
10800

Customer Tag No.:

Certificate Details

Number:	3068/20	Date of Issue:	22-Jul-2020	Expiry date:	22-Jul-2022
Material Details					
Production Order:	90161003	Material Code:	640300-SK-44	Cylinder No.:	A00692SK
Gas content:	5.52 M ³	Filling pressure:	145.0 bar	Valve:	CGA 660 SS
Cylinder Owner:	LINDE	Cylinder Material:	Spectra seal	Cylinder Size:	40 L

Laboratory Report

Analytical Result

Component	Normal Concentration	Analysis Result ¹	Uncertainty ²	Method of Analysis ³	Assay Date
Nitric Oxide	40.0 ppm	40.8 ppm	± 1% relative	(6) I-PB-352	15-Jul & 22-Jul-20
Other NOx impurity In Nitrogen		Less than 2.0 ppm			

Reference Standard used in Assay

Reference Standard	Cylinder number	Concentration	Expiry date:
Nitric Oxide In Nitrogen	7658395G	51.22 ± 0.41 ppm	30-Oct-2020

Analytical Instruments used in Assay

Instrument/Make/Model	Analytical Principle	Last Multipoint Calibration
FTIR Spectrometers Nicolet iSS0	FTIR-NO	19-Jun & 20-Jul-20

Recommend usage condition

Minimum utilization: 5% of actual content or before expire date whichever comes first.

Storage condition: Keep in well ventilation and secure area.

Comments

When reordering, please quote the material number

Note:

1. All results expressed in this report are on mole/mole basis, unless otherwise specified. The Assay of this Standard has been performed in accordance with the EPA Traceability Protocol EPA-600/R-12/531 for the Assay and Certification of Gaseous Calibration Standards using procedure G1
2. The reported expanded uncertainty is based on a standard uncertainty multiplied by a coverage factor k=2, providing a level of confidence of approximately 95%. The measurement of this material is traceable to the SI through the reference gas standard which is traceable to Swiss National Standard of Mass or other recognised national metrology institutes.
3. (1) Gas Chromatography, (2) Paramagnetic Oxygen Analyzer, (3) Electrochemical Oxygen Analyzer, (4) Electrochemical Moisture Analyzer, (5) Total Hydrocarbon Analyzer, (6) Other - Specified

Sukanya Parinyasoontorn

Signatory for and on behalf of Linde (Thailand) Co., Ltd.

Page 1 of 1

This report shall not be reproduced except in full

บริษัท ลินด์ (ประเทศไทย) จำกัด (มหาชน)

เลขที่ใบรับแจ้งการขึ้นทะเบียน 0107537000785

ชั้น 15 อาคารทาวเวอร์ เอ 2/3 หมู่ 14 ถนนบางนา-ตราด กม. 6.5 คลองแก้ว

อ.บางพลี จ.สมุทรปราการ 10540 โทรศัพท์ (66) 2338-6100 โทรสาร (66) 2338-6333

โรงงานเวลโกรว์: 105 หมู่ 5 คลองสกลนคร อ.บางพลีใหญ่ จ.สมุทรปราการ 24180

โทรศัพท์ (66) 38.570-479-93

โทรสาร (66) 38.570-323

Linde (Thailand) Public Company Limited 381/1, 01 October 2019

PLC, Registration no. 0107537000785

15th Floor, Bangna Tower A, 2/3 Moo 14, Bangna Trad KM. 6.5 Road, Bangkaew

Bangplee, Samutprakarn 10540, Tel (66) 2338-6100 Fax (66) 2338-6333

Wellgrow Plant: 105 Moo 5, T.Bangsamak, A.Bangpakong, Chachoengsao 24180

Thailand, Tel (66) 38.570-479-93

Fax (66) 38.570-323

THE LINDE GROUP

Linde

Certificate Of Analysis
Special Gases Mixture

Customer Details

Name:

Secot Co., Ltd.

Address:

239, Rimklongprapa Rd., Bangsue, Bangkok
10800

Customer Tag No.:

Certificate Details

Number:	3073/20	Date of Issue:	23-Jul-2020	Expiry date:	23-Jul-2022
Material Details					
Production Order:	90161004	Material Code:	614500-SK-44	Cylinder No.:	A00980SK
Gas content:	5.52 M ³	Filling pressure:	145.0 bar	Valve:	CGA 660 SS
Cylinder Owner:	LINDE	Cylinder Material:	Spectra seal	Cylinder Size:	40 L

Laboratory Report

Analytical Result

Component	Nominal Concentration	Analysis Result ¹	Uncertainty ²	Method of Analysis ³	Assay Date
Nitric Oxide	80.0 ppm	82.4 ppm	± 1% relative	(6) I-PB-352	16-Jul & 23-Jul-20
Other NOx impurity In Nitrogen		Less than 4.1 ppm			

Reference Standard
Nitric Oxide
In NitrogenReference Standard used in Assay
Cylinder number
2658395GConcentration
51.22 ± 0.41 ppmExpiry date:
30-Oct-2020Instrument/Make/Model
FTIR Spectrometers Nicolet iS50Analytical Instruments used in Assay
Analytical Principle
FTIR-NOLast Multipoint Calibration
19-Jun & 20-Jul-20

Recommend usage condition

Minimum utilization: 5% of actual content or before expire date whichever comes first.
Storage condition: Keep in well ventilation and secure area.

Comments

When reordering, please quote the material number

Note:

1. All results expressed in this report are on mole/mole basis, unless otherwise specified. The Assay of this standard has been performed in accordance with the EPA Traceability Protocol EPA-600/R-12/531 for the Assay and Certification of Gaseous Calibration Standards using procedure G1.
2. The reported expanded uncertainty is based on a standard uncertainty multiplied by a coverage factor k=2, providing a level of confidence of approximately 95%. The measurement of this material is traceable to the SI through the reference gas standard which is traceable to Swiss National Standard of Mass or other recognised national metrology institutes.
3. (1) Gas Chromatography, (2) Paramagnetic Oxygen Analyzer, (3) Electrochemical Oxygen Analyzer, (4) Electrochemical Moisture Analyzer, (5) Total Hydrocarbon Analyzer, (6) Other - Specified

Page 1 of 1

This report shall not be reproduced except in full

บริษัท ลินด์ (ประเทศไทย) จำกัด (มหาชน)

เลขที่ใบอนุญาตประกอบวิชาชีพ 0107537000185

ชั้น 15 อาคารทาวเวอร์ เอ 2/3 หมู่ 14 ถนนพหลโยธิน-ตราด กม. 6.5 ต.บางพลี

อ.บางพลี จ.สมุทรปราการ 10540 โทรศัพท์ (66) 2338-6100 โทรสาร (66) 2338-6333

โรงงานเวลโกรว์: 105 หมู่ 5 ต.บางพลี อ.บางพลี จ.สมุทรปราการ 24180

โทรศัพท์ (66) 38.570-479-93

โทรสาร (66) 38.570-323

Sukanya Parinyasoonorn

Signatory for and on behalf of Linde (Thailand) Co., Ltd.

PB-002/1006

Linde (Thailand) Public Company Limited 1/1, 01 October 2019

PIC Registration no 0107537000185

15th Floor, Bangna Tower A, 2/3 Moo 14, Bangna Trad KM. 6.5 Road, Bangkaew

Bangplee, Samutprakarn 10540, Tel (66) 2338-6100 Fax (66) 2338-6333

Wellgrow Plant, 105 Moo 5, 1, Bangsamak, A, Bangpakong, Chachoengsao 24180

Thailand, Tel (66) 38.570-479-93

Fax (66) 38.570-323



Request Service No. 076/63

Page 1 of 3

Calibration Certificate

Nomenclature : Brand : Mettler Toledo Type : Top-Loading Electronic Balance

Model : AG 245 Serial No. : 1117293916 (198129-0)

Submitted by : Laboratory of SECOT CO., LTD.

Location of Calibration : BAL Room , 6th Floor, Secot Co., Ltd.

Calibration range : 0 – 200 g Scale division : 0.00001 g (41 g)/0.0001 g (210 g)

Calibration date : May 28, 2020

Reference Standard No. M1912055S, M200162, M190247

Traceable to : Metrological Center SCI ECO Services Co., Ltd, THAI CALIBRATION SERVICES Co., Ltd.

Ambient Condition : Temperature 24.92-25.42 °C

Humidity 46.1 - 51.0 % RH

Calibrated By :

Suratwadee Chaiyatham

(Miss Suratwadee Chaiyatham)

Testing Officer

Date : 29 / 05 / 2020

Approved By :

Narisa Poowasanpetk

(Miss Narisa Poowasanpetk)

Chief of Technical Management

Date : 29 / 05 / 2020

Issued Date : May 29, 2020

Measurement Report

Request Service No. 076/63

Page 2 of 3

Description : Brand : Mettler Toledo

Type : Top-Loading Electronic Balance

Model : AG 245

Serial No. : 1117293916 (198129-0)

Calibration range : 0 – 200 g

Scale division : 0.00001 g (41 g)/0.0001 g (210 g)

Calibration date : May 28,2020

Ambient Condition : Temperature 24.92-25.42 °C Relative humidity 46.1-51.0 % RH

Measurement data :

1. Repeatability of Reading :

Load (g)	Standard Deviation of Reading (g)	Maximum Difference between Successive Reading (g)
50	0.00005	0.0001
100	0.00005	0.0001
150	0.00008	0.0002
200	0.00017	0.0006

2. Off-Center Loading :

A Mass of 50.0000 g was placed and moved to various position on the pan.

Unit : g

Center	Front	Left	Back	Right	Center	Maximum Difference
50.00054	50.00076	50.00060	50.00022	50.00030	50.00054	0.00032

Issued Date : May 29,2020

3. Departure from Nominal Value :

Reading (g)	Correction (g)	Uncertainty (+/- g)
0	0.000000	± 0.000033
0.5	+0.000003	± 0.000034
1	+0.000029	± 0.000024
10	-0.000130	± 0.000035
20	-0.000222	± 0.000050
40	-0.000471	± 0.000077
60	-0.00060	± 0.00016
80	-0.00072	± 0.00014
100	-0.00097	± 0.00017
120	-0.00100	± 0.00019
140	-0.00110	± 0.00021
160	-0.00121	± 0.00024
180	-0.00135	± 0.00026
200	-0.00164	± 0.00028

Calibrated by :

Suratwadee Chaiyatham

Approved By :

Narisa Poowasanpetch

(Miss Suratwadee Chaiyatham)

(Miss Narisa Poowasanpetch)

Testing Officer

Chief of Technical Management

Date : 29 / 05 / 2020

Date : 29 / 05 / 2020

Issued Date : May 29, 2020



Request Service No.074/63

Page 1 of 3

Calibration Certificate

Nomenclature : Brand : Mettler Toledo Type : Top-Loading Electronic Balance

Model : BSA224S-CW Serial No. : 32191636

Submitted by : Laboratory of SECOT CO., LTD.

Location of Calibration : BAL Room , 6th Floor, Secot Co., Ltd.

Calibration range : 0 – 200 g

Scale division : 0.0001 g (220 g)

Calibration date : May 25,2020

Reference Standard No. M200162, M190247

Traceable to : Metrological Center SCI ECO Services Co.,Ltd.

Ambient Condition : Temperature 25.56-26.33 °C

Humidity 48.8-51.1 % RH

Calibrated By : *Khemchuda Insorn*

Approved By : *Narisa Poowasanpetch*

(Miss Khemchuda Insorn)

(Miss Narisa Poowasanpetch)

Testing Officer

Chief of Technical Management

Date : *26/05/2020*

Date : *26/05/2020*

Issued Date : May 26,2020

Measurement Report

Request Service No.074/63

Page 2 of 3

Description : Brand : Sartorius

Type : Top-Loading Electronic Balance

Model : BSA224S-CW

Serial No. : 32191636

Calibration range : 0 – 200 g

Scale division : 0.0001 g (220 g)

Calibration date : May 25,2020

Ambient Condition : Temperature 25.56-26.33 °C Relative humidity 48.8-51.1 % RH

Measurement data :

1. Repeatability of Reading :

Load (g)	Standard Deviation of Reading (g)	Maximum Difference between Successive Reading (g)
50	0.00006	0.0002
100	0.00007	0.0002
150	0.00004	0.0001
200	0.00005	0.0001

2. Off-Center Loading :

A Mass of 50.0000 g was placed and moved to various position on the pan.

Unit : g

Center	Front	Left	Back	Right	Center	Maximum Difference
49.99976	49.99978	49.99974	49.99986	49.99982	49.99986	0.00012

Issued Date : May 26,2020

3. Departure from Nominal Value :

Reading (g)	Correction (g)	Uncertainty (+/- g)
0	0.00000	± 0.00008
1	- 0.00002	± 0.00008
5	+ 0.00009	± 0.00011
10	+ 0.00035	± 0.00008
20	+ 0.00009	± 0.00010
40	+0.00013	± 0.00010
60	+0.00020	± 0.00012
80	+0.00020	± 0.00014
100	+0.00015	± 0.00016
120	+0.00014	± 0.00018
140	+0.00022	± 0.00020
160	+0.00009	± 0.00022
180	+0.00010	± 0.00025
200	+0.00006	± 0.00027

Calibrated by : Khemchuda InsornApproved By : Narisa Poowasanpetch

(Miss Khemchuda Insorn)

(Miss Narisa Poowasanpetch)

Testing Officer

Chief of Technical Management

Date : 26/05/2020Date : 26/05/2020

Issued Date : May 26, 2020



Bangkok High Lab Co.,Ltd.

4/176 Soi Ladplakao 66, Ladplakao Rd., Anusawari, Bangkok, Bangkok 10220

Tel: (662) 971-5800

Website: www.bangkokhighlab.com

Fax: (662) 971-5300

E-mail: info@bangkokhighlab.com



NSC-TISI-TIS 17025
CALIBRATION 0366

CERTIFICATE OF CALIBRATION

Certificate No : S2020/187

Page : 1/6

Order : 460/2020

Customer : Secot Co., Ltd.
Address : 239 Rimklongprapa Rd., Bangsue, Bangkok 10800
Instrument : UV/VIS spectrophotometer
Manufacture : Thermo Scientific
Model : Genesys 10s uv-vis
Serial Number : 2L9N349007
Environment : Temperature (25.2 - 24.1) °C
: Humidity (56 - 59) %RH
Received Date : November 3, 2020
Calibration Date : November 3, 2020
Issue Date : November 11, 2020
Calibrate Status : No Adjustment
Calibration Area : Customer area
Roomname : Laboratory Room of Secot Co., Ltd.
Calibrated By : Wisamuchai Ngiruned

(Mr. Wisanuchai Ngiruned)

Calibration Engineer

Approved By

: Wanchai Meesiri

(Mr. Wanchai Meesiri)

Manager



Bangkok High Lab Co.,Ltd.

4/176 Soi Ladplakao 66, Ladplakao Rd., Anusavari, Bangkok, Bangkok 10220

Tel: (662) 971-5800

Website: www.bangkokhighlab.com

Fax: (662) 971-5300

E-mail: info@bangkokhighlab.com



NSC-TISI-TIS 17025
CALIBRATION 0366

Certificate No : S2020/187

Page : 2/6

Order No : 460/2020

1. Photometric Accuracy

CRMs: Neutral Density Glass Filters

CRMs Serial Number: 10563

Traceability: Traceable to NIST, U.S.A. through Neutral density filters NIST SRM 930e & 1930, Double Aperture method through Starna certificate report 77969

Spectral slit width : 1.80 nm

1.1 Reading scale at 420.0 nm

Filter STDs (Abs) Certificate	Average Measured Value (A)	Correction (A)	Uncertainty \pm (A)
0.0000	0.000	0.0000	0.0028
0.5571	0.557	0.0001	0.0048
1.0732	1.070	0.0033	0.0068
2.2232	2.228	-0.0043	0.0110

1.2 Reading scale at 440.0 nm

Filter STDs (Abs) Certificate	Average Measured Value (A)	Correction (A)	Uncertainty \pm (A)
0.0000	0.000	0.0000	0.0028
0.5470	0.548	-0.0010	0.0048
1.0465	1.046	0.0005	0.0046
2.1566	2.160	-0.0031	0.0064

1.3 Reading scale at 465.0nm

Filter STDs (Abs) Certificate	Average Measured Value (A)	Correction (A)	Uncertainty \pm (A)
0.0000	0.000	0.0000	0.0028
0.4969	0.497	-0.0002	0.0040
0.9619	0.964	-0.0021	0.0040
2.0055	2.008	-0.0020	0.0060



Bangkok High Lab Co., Ltd.

4/176 Soi Ladplakao 66, Ladplakao Rd., Anusawari, Bangkok, Bangkok 10220

Tel: (662) 971-5800

Website: www.bangkokhighlab.com

Fax: (662) 971-5300

E-mail: info@bangkokhighlab.com



NSC-TISI-TIS 17025
CALIBRATION 0366

Certificate No : S2020/187

Page : 3/6

Order No : 460/2020

1.4 Reading scale at 546.1nm

Filter STDs (Abs) Certificate	Average Measured	Correction (A)	Uncertainty \pm (A)
0.0000	0.000	0.0000	0.0028
0.5109	0.511	-0.0001	0.0029
0.9748	0.971	0.0037	0.0030
2.0254	2.028	-0.0027	0.0060

1.5 Reading scale at 590 nm

Filter STDs (Abs) Certificate	Average Measured Value (A)	Correction (A)	Uncertainty \pm (A)
0.0000	0.000	0.0000	0.0028
0.5391	0.540	-0.0009	0.0034
1.0157	1.012	0.0038	0.0048
2.0647	2.070	-0.0051	0.0067

1.6 Reading scale at 635 nm

Filter STDs (Abs) Certificate	Average Measured Value (A)	Correction (A)	Uncertainty \pm (A)
0.0000	0.000	0.0000	0.0028
0.5234	0.528	-0.0046	0.0058
0.9683	0.967	0.0012	0.0060
1.9526	1.955	-0.0024	0.0068



Bangkok High Lab Co.,Ltd.

4/176 Soi Ladplakao 66, Ladplakao Rd., Anusawari, Bangkhen, Bangkok 10220

Tel: (662) 971-5800

Website: www.bangkokhighlab.com

Fax: (662) 971-5300

E-mail: info@bangkokhighlab.com



NSC-TISI-TIS 17025
CALIBRATION 0366

Certificate No : S2020/187

Page : 4/6

Order No : 460/2020

2. Photometric Accuracy

CRMs: Potassium Dichromate in Perchloric acid

CRMs Serial Number: 16792

Blank Serial Number: 16708

Traceability: Traceable to NIST, U.S.A. through crystalline potassium dichromate NIST SRM 935a through Starna certificate report 77966

Spectral slit width : 1.80 nm

Wave Length	Certificate (Abs)	Average Measured Value (A)	Correction (A)	Uncertainty \pm (A)
235	0.0000	0.000	0.0000	0.0050
	0.7295	0.729	0.0005	0.0056
257	0.0000	0.000	0.0000	0.0050
	0.8499	0.850	0.0000	0.0054
313	0.0000	0.000	0.0000	0.0050
	0.2853	0.288	-0.0027	0.0054
350	0.0000	0.000	0.0000	0.0050
	0.6316	0.632	-0.0004	0.0054

3. Wavelength Accuracy

Spectral slit width : 1.80 nm

3.1 CRMs: Holmium Glass Filter

CRMs Serial Number: 10763

Traceability: Traceable to NIST Holmium oxide filter NIST SRM 2034, through Starna certificate report 77965

Filter STDs (nm) Certificate	Average Measured Value (nm)	Correction (nm)	Uncertainty \pm (nm)
241.74	240.6	1.14	0.12
279.44	279.0	0.44	0.12
287.98	287.8	0.18	0.12
334.10	334.0	0.10	0.12
361.00	360.8	0.20	0.12
418.61	418.8	-0.19	0.12
453.63	453.8	-0.17	0.12
460.05	460.4	-0.35	0.12
536.66	536.6	0.06	0.12
637.98	638.0	-0.02	0.12



Bangkok High Lab Co.,Ltd.

4/176 Soi Ladplakao 66, Ladplakao Rd., Anusawari, Bangkok, Bangkok 10220

Tel: (662) 971-5800

Fax: (662) 971-5300

Website: www.bangkokhighlab.com

E-mail: info@bangkokhighlab.com



NSC-TISI-TIS 17025
CALIBRATION 0366

Certificate No : S2020/187

Page : 5/6

Order No : 460/2020

3.2 CRMs: Didymium Glass Filter

CRMs Serial Number: 10764

Traceability: Traceable to NIST Didymium filter NIST SRM 2034, through Starna certificate report 77964

Filter STDs (nm) Certificate	Average Measured Value (nm)	Correction (nm)	Uncertainty ± (nm)
585.29	585.6	-0.31	0.12
684.49	684.8	-0.31	0.12
740.18	740.6	-0.42	0.12
748.48	749.0	-0.52	0.12
807.03	807.6	-0.57	0.12
879.27	879.8	-0.53	0.12

4. *Stray Light

CRMs: Potassium Chloride aqueous solution

CRMs Serial Number: 14912

Blank Serial Number: 14958

Traceability: Traceable to NIST, U.S.A. crystalline potassium chloride NIST SRM2032, through Starna certificate report 77967

Spectral slit width : 1.80 nm

Wavelength (nm)	Certificate	Average Measured
201.07	>2A	#N/A
201.07	<1%T	#N/A



Bangkok High Lab Co.,Ltd.

4/176 Soi Ladplakao 66, Ladplakao Rd., Anusawari, Bangkok, Bangkok 10220

Tel: (662) 971-5800

Fax: (662) 971-5300

Website: www.bangkokhighlab.com

E-mail: info@bangkokhighlab.com



NSC-TISI-TIS 17025
CALIBRATION 0366

Certificate No : S2020/187

Page : 6/6

Order No : 460/2020

5.*Spectral Resolution

CRMs: Toluene in Hexane

CRMs Serial Number: 14812

Blank Serial Number: 14803

Traceability: Traceable to toluene in hexane NIST SRM2034, through Starna certificate report 77968

Spectral slit width (nm)	Abs Ratio
0.5	#N/A
1.0	#N/A
1.5	#N/A
2.0	#N/A
3.0	#N/A

Note : * "Not TISI Accredited" in this certificate have been included for completeness

Remark:

1. Calibrate Method
 - 1.1 Photometric accuracy: In-house method W-SER-001 based on ASTM E925-02 and ASTM E275-01
 - 1.2 Wavelength accuracy: In-house method W-SER-001 based on ASTM E925-02 and ASTM E275-01
 - 1.3 Stray light: Measuring the CRMs in both absorbance and transmittance unit at wavelength 201.23 nm. Base on European Pharmacopoeia V.6.19.3 1984
 - 1.4 Spectral resolution: Measuring the CRMs. The maximum absorbance values were read at closest to 268.7nm and the minimum absorbance values were read at closest 267.0nm. Refer to European Pharmacopoeia V.6.19.3 1984
2. N/A = not available.
3. Uncertainty of Measurement: The reported uncertainty of measurement was based on standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%
4. This result of calibration was found accurate as shown on date and place of calibration only.
5. This report will certify of calibrated equipment only.

- End of Report -



National Food Institute, Ministry of Industry, Thailand

2008 Soi Arun Ammarin 36, Arun Ammarin Rd., Bangyeekhan, Bangphlad, Bangkok 10700 Thailand
Tel : +66 (0) 2422 8688 Fax : +66 (0) 2422 8558 Website : www.nfi.or.th E-mail : cal@nfi.or.th



Calibration Certificate

Certificate No.: 2004097-001-01
Client name: SECOT CO., LTD.
Address: 239 Rimklongprapa Road, Bangsue,
Bangsue, Bangkok 10800

Page 1 of 3

Equipment: CHAMBER (Hot Air Oven)

Manufacturer: BINDER

Model: ED 53

Serial No.: 01-27152

ID No.: N/A

Order No.: 2004097

Operation No.: 2004097-001

Date of Receipt: 3 August 2020

Date of Calibration: 3 August 2020

Calibrated by Mr.Pheraphat Tuanjit
Scientist

Approved by

(Mr.Manas Somsak)

Senior Analyst, Division of Calibration Laboratory

Date of Issue: 5 August 2020

Responsible for the Technical Management Team

The uncertainties are for a confidence probability of approximately 95 %.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

Calibration Report

Certificate No.: 2004097-001-01
Equipment: CHAMBER (Hot Air Oven)
Model: ED 53 Serial No.: 01-27152
Resolution: 1 °C ID No.: N/A
Manufacturer: BINDER

Date of Calibration: 3 August 2020

Page 2 of 3

Location: Walkway, SECOT CO., LTD.
Environment Condition: Ambient Temperature (32 ± 1) °C
Relative Humidity (55 ± 5) %
Line Voltage (229 ± 1) Volt

Condition of this results of Calibration:

- This instrument was calibrated by insert 9 standard thermometer into its chamber and calibration according to W-TE-014 Based on TLAS G-20-1/02-08 (E): Guidelines for Calibration and Checks of Temperature Controlled Enclosures.
- The temperature scale used was based on ITS - 90.
- All data show below were final values and the initial data may be obtained upon request.
- Reference Standard Instrument :

Instrument	Model	Serial No./ID No.	Certificate No.	Due Date	Through
Digital Thermometer with sensor	34972A	MY49020132	TE 630168-01	27 January 2021	NATIONAL FOOD INSTITUTE
	RTD	CH#301-309/ RTD#301-309			

- This certificate is traceable to International System of Units (SI Units).
- This certificate was certified only for the instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.
- Condition of Calibrated item : Good

UUC Description :

Time of Record 1 Hour 9 Minute At 104 °C
Fresh air Damper ☐ Open Position ☐
☒ Close
☐ Not Available

7. Result of Calibration : ☒ Without adjustment ☐ After adjustment

an

Calibration Report

Certificate No.: 2004097-001-01
Equipment: CHAMBER (Hot Air Oven)
Model: ED 53 **Serial No.:** 01-27152
Resolution: 1 °C **ID No.:** N/A
Manufacturer: BINDER

Date of Calibration: 3 August 2020

Page 3 of 3

Calibration point: 104, 110 and 180 °C

Calibration result:

Calibration Condition	Temperature (°C)	Relative Humidity (%)	Line Voltage (Volt)
MIN	31.3	50	228.8
MAX	32.1	60	230.0

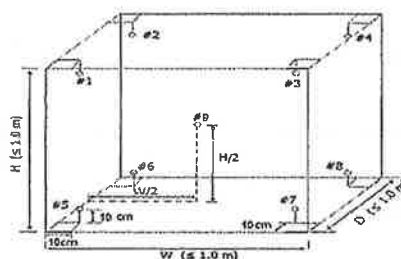


Table1 : Reporting of Temperature

Calibration point (°C)	Measured Temperature (°C) @ Sensor No. (Sensor No.9 is REF)									Uncertainty ± (°C)
	# 1	# 2	# 3	# 4	# 5	# 6	# 7	# 8	# 9	
104	104.55	104.43	104.69	104.23	103.83	103.04	103.63	103.53	103.20	0.81
110	110.64	110.59	110.76	110.29	109.91	109.07	109.70	109.58	109.24	0.82
180	180.15	180.46	179.88	179.84	179.81	178.73	179.45	179.73	178.97	0.91

Table 2 : Reporting of Characterization Result

UUC* Setting (°C)	UUC* reading (°C)			Stability ± (°C)	Uniformity (°C)	Overall Variation (°C)
	MIN	MAX	Average			
104	104	104	104	0.20	1.62	1.96
110	110	110	110	0.22	1.63	1.95
177	177	177	177	0.31	1.75	2.28

Note The quoted uncertainty include " Stability " and " Loading effect (20% of Temp Uniformity) "

UUC* = Unit Under Calibration

Stability = One-half of the greatest maximum difference of measured temperatures at any one sensors, for at least half an hour after reaching steady state.

Uniformity = The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time.

Overall Variation = The difference of the maximum and minimum measured temperatures throughout observation time.

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor $k=2$, providing a level of confidence of approximately 95 %.

----- End -----

Signature



National Food Institute, Ministry of Industry, Thailand

2008 Sol 36, Arun Amarin Road, Bang Yi Khan Subdistrict, Bang Phlat District, Bangkok 10700, Thailand.

Tel : +66 (0) 2422 8688 Fax : +66 (0) 2422 8545 Website : www.nfi.or.th E-mail : cal@nfi.or.th



Calibration Certificate

Certificate No.: 2103125-001-01
Client name: SECOT CO.,LTD.
Address: 239 Rimklongprapa Road, Bangsue,
Bangsue, Bangkok 10800

Page 1 of 3

Equipment: CHAMBER (Hot Air Oven)

Manufacturer: MEMMERT

Model: UF 55

Serial No.: B213.0295

ID No.: N/A


Order No.: 2103125

Operation No.: 2103125-001

Date of Receipt: 4 June 2021

Date of Calibration: 4 June 2021

Calibrated by Mr.Pheraphat Tuanjit
Scientist

Approved by 
(Mr.Nuttapol Niyomchat)
Expert, Division of Calibration Laboratory

Date of Issue: 10 June 2021

Responsible for the Technical Management Team

The uncertainties are for a confidence probability of approximately 95 %.

This Certificate is Issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.



National Food Institute, Ministry of Industry, Thailand

2008 Soi 36, Arun Amarin Road, Bang Yi Khan Subdistrict, Bang Phlat District, Bangkok 10700, Thailand.
Tel : +66 (0) 2422 8688 Fax : +66 (0) 2422 8545 Website : www.nfi.or.th E-mail : cal@nfi.or.th



Calibration Report

Certificate No.: 2103125-001-01
Equipment: CHAMBER (Hot Air Oven)
Model: UF 55 Serial No.: B213.0295
Resolution: 0.1 °C ID No.: N/A
Manufacturer: MEMMERT

Date of Calibration: 4 June 2021

Page 2 of 3

Location: Walkway Laboratory, SECOT CO.,LTD.
Environment Condition: Ambient Temperature (35 ± 1) °C
Relative Humidity (50 ± 6) %
Line Voltage (217 ± 3) Volt

Condition of this results of Calibration:

- This instrument was calibrated by insert 9 standard thermometer into its chamber and calibration according to W-TE-014 Based on TLAS G-20-1/02-08 (E): Guidelines for Calibration and Checks of Temperature Controlled Enclosures.
 - The temperature scale used was based on ITS - 90.
 - All data show below were final values and the Initial data may be obtained upon request.

2. Reference Standard Instrument :

Instrument	Model	Serial No./ID No.	Certificate No.	Due Date	Through
Digital Thermometer with sensor	34972A	MY49010812	TE 640401-01	28 April 2022	NATIONAL FOOD INSTITUTE
	RTD	CH#101-109/ RTD#101-109			

- This certificate is traceable to International System of Units (SI Units).
- This certificate was certified only for the Instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.
- Condition of Calibrated item : Good

UUC Description :

Time of Record 1 Hour 9 Minute At 80.0 °C
Fresh air Damper ☐ Open Position ☐
☒ Close
☐ Not Available

7. Result of Calibration : ☒ Without adjustment ☐ After adjustment

P. Jaenghanvit
10 June 2021

Calibration Report

Certificate No.: 2103125-001-01
Equipment: CHAMBER (Hot Air Oven)
Model: UF 55 **Serial No.:** B213.0295
Resolution: 0.1 °C **ID No.:** N/A
Manufacturer: MEMMERT

Date of Calibration: 4 June 2021

Page 3 of 3

Calibration point: 80.0, 104.0 and 180.0 °C

Calibration result:

Calibration Condition	Temperature (°C)	Relative Humidity (%)	Line Voltage (Volt)
MIN	33.7	44	214.0
MAX	35.8	55	220.0

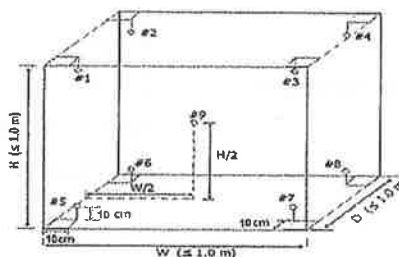


Table1 : Reporting of Temperature

Calibration point (°C)	Measured Temperature (°C) @ Sensor No. (Sensor No.9 is REF)									Uncertainty ± (°C)
	# 1	# 2	# 3	# 4	# 5	# 6	# 7	# 8	# 9	
80.0	79.85	79.78	80.07	79.89	79.80	79.75	79.59	80.13	79.85	0.46
104.0	104.03	103.90	104.29	104.08	103.99	103.91	103.63	104.36	104.05	0.53
180.0	180.10	179.71	180.35	180.02	180.06	179.82	179.42	180.91	179.99	0.90

Table 2 : Reporting of Characterization Result

UUC* Setting (°C)	UUC* reading (°C)			Stability ± (°C)	Uniformity (°C)	Overall Variation (°C)
	MIN	MAX	Average			
80.0	80.0	80.0	80.0	0.06	0.28	0.68
104.0	104.0	104.0	104.0	0.10	0.43	0.92
180.0	180.0	180.0	180.0	0.14	0.92	1.75

Note The quoted uncertainty include " Stability " and " Loading effect (20% of Temp Uniformity) "

UUC* = Unit Under Calibration

Stability = One-half of the greatest maximum difference of measured temperatures at any one sensors, for at least half an hour after reaching steady state.

Uniformity = The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time.

Overall Variation = The difference of the maximum and minimum measured temperatures throughout observation time.

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k= 2, providing a level of confidence of approximately 95 %.

----- End -----

P. Jongsakulrat
10 June 2021



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG, BANGKOK 10250
TEL. 0-2717-3000-24 FAX. 0-2719-9484



Certificate of Calibration

Certificate No. : 20H2194

Page : 1 of 2

Equipment : Digital Thermo-Hygrometer
Manufacturer: Digicon
Model : HT-776
Serial No.: Q8588B6
ID No.: 214400

This certificate may not be reproduced other than in full,
except with the prior written approval of the head of
Corporate Services 3: Equipment Calibration and Testing Services.

Condition As-Received: Used Item

Received Date: 08 September 2020

Calibration Date: 10 September 2020
to 17 September 2020

Reference: 2009-0153DN

Submitted by: Secot Co.,Ltd.

Ambient Temperature: $(25 \pm 3) ^\circ\text{C}$

Relative Humidity: $(50 \pm 20) \%$

239 Rimklongprapa Road., Bangsue, Bangkok 10800

Procedure used: Calibration were conducted using in-house calibration procedure CP-H03 according to comparison with standard chilled mirror sensor for humidity measurement function and comparison with standard temperature probe for temperature measurement function into humidity / temperature chamber.

Condition of this result of calibration

1. Reference standards instruments :

<u>Instrument</u>	<u>Model</u>	<u>Serial No.</u>	<u>Certificate No.</u>	<u>Due Date</u>
1) Chilled Mirror Hygrometer	Dew Master	44730	18203	07 Apr 2021
2) Handheld Thermometer With Sensor	1523	3240076	201203	13 Feb 2021

2. The certificate is valid only to the item calibrated on date and place of calibration.

3. This Certification is traceable to the International System of Unit maintained at:-

- National Institute of Standards and Technology (NIST) , The United States of America
- National Institute of Metrology Thailand (NIMT)

Calibrated by : Kralpop Onrat
Issue Date : 18 September 2020

Approved Signatory :

[✓] Chakrit Waewanjua

[] Ponthippa Tameyakul

[] Pitak Srimongkol

B 0242537



Cert. No.: 20H2194

Page.: 2 of 2

This instrument was connected with humidity/temperature probe Serial No. Q858886.

Result of Calibration:-

Without Adjustment

Function:

Humidity measurement.

Reference Temperature (°C)	Standard Humidity (%R.H.)	UUC* Reading (%R.H.)	Error (%R.H.)	Uncertainty of Measurement (±%R.H.)
25.0	40.1	42.28	2.18	1.3
25.0	50.1	50.82	0.72	1.6
25.0	60.0	59.17	-0.83	1.6

Result of Calibration:-

Without Adjustment

Function:

Temperature measurement.

Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty of Measurement (±°C)
20.030	20.20	0.170	0.42
25.027	25.25	0.223	0.42
30.004	30.28	0.276	0.42

UUC* : Unit Under Calibration

The reported uncertainty of measurement was base on standard uncertainty multiplied by coverage factor $k = 2.00$, providing confidence level approximately 95%.

-o0o-

a 1019823

CAL

Calibratech Co.,Ltd.

7/106-7 Moo 2, Sukhaphrachasan 3 Rd., Bangpood, Pakkred, Nonthaburi 11120

Tel.(02) 964-6211 Fax.(02) 964-5155, e-mail : calibratech-cal@yahoo.com, calibratech-cn@hotmail.com



NSC-TISI-TIS 17025
CALIBRATION 0030

Certificate of Calibration

Certificate No. : 64-420032-1

Page : 1 of 2

Submitted by : Secot Co.,Ltd.

239 RimKlongprapa Road, Bangsue, Bangkok 10800 Thailand

Equipment : pH Meter with electrode

pH meter

Manufacturer : Mettler Toledo

Model : Seven2Go S2

Range : N/A pH

Resolution : 0.01 pH

Serial No. : B924795409

ID No. : N/A

Electrode

Model : InLab Expert Go-ISM

Serial No. : 0436432

Environment : Ambient Temperature : $(25 \pm 2) ^\circ \text{C}$

Relative Humidity : $(50 \pm 15) \%$

Date of Received : 15 February 2021

Date of Calibration : 17 February 2021

Date of Issue : 17 February 2021

Calibrated by : Bunjerd Masri

Calibration Method : In-house method CAL-M4201 direct measurement by using standard voltage calibrator and using certified reference material (CRM)


Reference Standard Instruments : This certification is traceable to the International System of Units

1. Multiproduct Calibrator

<u>ID No.</u>	<u>Cert. No.</u>	<u>Due Date</u>	<u>Traceability</u>
400005	E1U190739	31 Aug 2021	National Institute of Metrology Thailand (NIMT)

2. Standard Buffer Solution

<u>pH</u>	<u>Cert. No.</u>	<u>Lot No.</u>	<u>Exp. Date</u>	<u>Traceability</u>
4.004	61208711	684575	10 Apr 2021	CPA chem
6.985	61191143	684576	10 Apr 2021	CPA chem
9.963	61208865	684577	10 Apr 2021	CPA chem

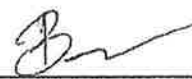Approved by : 
(Bunjerd Masri)
Supervisor

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full except with the prior written approval of the Calibratech Co.,Ltd.



CAL

Calibratech Co.,Ltd.

7/106-7 Moo 2, Sukhprachasan 3 Rd., Bangpood, Pakkred, Nonthaburi 11120

Tel.(02) 964-6211 Fax.(02) 964-5155, e-mail : calibratech.cal@yahoo.com, calibratech.cal@hotmail.com

Certificate of Calibration

Certificate No. : 64-420032-1

Page : 2 of 2

Result of Calibration :

UUC Condition As-Received : Good

Function : Electrical measurement

pH meter

Performing standard curve by Multiproduct Calibrator at pH (4,7,10)

Adjustment Curve at nominal pH	Applied Voltage (mV)	Nominal Value (pH)	UUC Reading		Correction (mV)	Uncertainty (± mV)
			(pH)	(mV)		
4, 7, 10	177.4800	4	4.00	177	0	0.58
	0.0000	7	7.00	0	0	0.58
	-177.4800	10	10.00	-178	1	0.58

Function : pH meter with electrode

Performing a three - buffer standard curve using buffer nominal pH (4,7,10)

Adjustment Curve at nominal pH	Standard Buffer (pH)	UUC Reading (pH)	Correction (pH)	Uncertainty (± pH)
4, 7, 10	4.004	4.00	0.00	0.011
	6.985	7.00	-0.01	0.020
	9.963	10.00	-0.04	0.053

Remark

UUC : Unit Under Calibration

This result of calibration was found accurate as shown on date and place of calibration only.

This reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%

- o0o -





Request Service No. 076/63

Page 1 of 3

Calibration Certificate

Nomenclature : Brand : Mettler Toledo Type : Top-Loading Electronic Balance

Model : AG 245 Serial No. : 1117293916 (198129-0)

Submitted by : Laboratory of SECOT CO., LTD.

Location of Calibration : BAL Room , 6th Floor, Secot Co., Ltd.

Calibration range : 0 – 200 g Scale division : 0.00001 g (41 g)/0.0001 g (210 g)

Calibration date : May 28, 2020

Reference Standard No. M1912055S, M200162, M190247

Traceable to : Metrological Center SCI ECO Services Co., Ltd, THAI CALIBRATION SERVICES Co., Ltd.

Ambient Condition : Temperature 24.92-25.42 °C

Humidity 46.1 - 51.0 % RH

Calibrated By :

Suratwadee Chaiyatham

(Miss Suratwadee Chaiyatham)

Testing Officer

Date : 29 / 05 / 2020

Approved By :

Narisa Poowasanpetch

(Miss Narisa Poowasanpetch)

Chief of Technical Management

Date : 29 / 05 / 2020

Issued Date : May 29, 2020

Measurement Report

Request Service No. 076/63

Page 2 of 3

Description : Brand : Mettler Toledo

Type : Top-Loading Electronic Balance

Model : AG 245

Serial No. : 1117293916 (198129-0)

Calibration range : 0 – 200 g

Scale division : 0.00001 g (41 g)/0.0001 g (210 g)

Calibration date : May 28,2020

Ambient Condition : Temperature 24.92-25.42 °C Relative humidity 46.1-51.0 % RH

Measurement data :

1. Repeatability of Reading :

Load (g)	Standard Deviation of Reading (g)	Maximum Difference between Successive Reading (g)
50	0.00005	0.0001
100	0.00005	0.0001
150	0.00008	0.0002
200	0.00017	0.0006

2. Off-Center Loading :

A Mass of 50.0000 g was placed and moved to various position on the pan.

Unit : g

Center	Front	Left	Back	Right	Center	Maximum Difference
50.00054	50.00076	50.00060	50.00022	50.00030	50.00054	0.00032

Issued Date : May 29,2020

3. Departure from Nominal Value :

Reading (g)	Correction (g)	Uncertainty (+/- g)
0	0.000000	± 0.000033
0.5	+0.000003	± 0.000034
1	+0.000029	± 0.000024
10	-0.000130	± 0.000035
20	-0.000222	± 0.000050
40	-0.000471	± 0.000077
60	-0.00060	± 0.00016
80	-0.00072	± 0.00014
100	-0.00097	± 0.00017
120	-0.00100	± 0.00019
140	-0.00110	± 0.00021
160	-0.00121	± 0.00024
180	-0.00135	± 0.00026
200	-0.00164	± 0.00028

Calibrated by :

Suratwadee Chaiyatham

Approved By

Narisa Poowasanpetch

(Miss Suratwadee Chaiyatham)

(Miss Narisa Poowasanpetch)

Testing Officer

Chief of Technical Management

Date : 29 / 05 / 2020

Date : 29 / 05 / 2020

Issued Date : May 29, 2020

CAL

Calibratech Co.,Ltd.

7/106-7 Moo 2, Sukhprachasan 3 Rd., Bangpood, Pakkred, Nonthaburi 11120

Tel.(02) 964-6211 Fax.(02) 964-5155, e-mail : calibratech.cal@yahoo.com, calibratech.enk@hotmail.com



NSC-TISI-TIS17025
CALIBRATION 0030

Certificate of Calibration

Certificate No. : 64-420032-1

Page : 1 of 2

Submitted by : Secot Co.,Ltd.

239 RimKlongprapa Road, Bangsue, Bangkok 10800 Thailand

Equipment : pH Meter with electrode

pH meter

Manufacturer : Mettler Toledo

Model : Seven2Go S2

Range : N/A

pH

Resolution : 0.01 pH

Serial No. : B924795409

ID No. : N/A

Electrode

Model : InLab Expert Go-ISM

Serial No. : 0436432

Environment :

Ambient Temperature : $(25 \pm 2) ^\circ \text{C}$

Relative Humidity : $(50 \pm 15) \%$

Date of Received : 15 February 2021

Date of Calibration : 17 February 2021

Date of Issue : 17 February 2021

Calibrated by : Bunjerd Masri

Calibration Method : In-house method CAL-M4201 direct measurement by using standard voltage calibrator and using certified reference material (CRM)

Reference Standard Instruments : This certification is traceable to the International System of Units

1. Multiproduct Calibrator

<u>ID No.</u>	<u>Cert. No.</u>	<u>Due Date</u>	<u>Traceability</u>
400005	E1U190739	31 Aug 2021	National Institute of Metrology Thailand (NIMT)

2. Standard Buffer Solution

<u>pH</u>	<u>Cert. No.</u>	<u>Lot No.</u>	<u>Exp. Date</u>	<u>Traceability</u>
4.004	61208711	684575	10 Apr 2021	CPA chem
6.985	61191143	684576	10 Apr 2021	CPA chem
9.963	61208865	684577	10 Apr 2021	CPA chem

Approved by :

(Bunjerd Masri)

Supervisor

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full except with the prior written approval of the Calibratech Co.,Ltd.

CAL-F0031-03



CAL

Calibratech Co.,Ltd.

7/106-7 Moo 2, Sukhprachasan 3 Rd., Bangpood, Pakkred, Nonthaburi 11120

Tel.(02) 964-6211 Fax.(02) 964-5155, e-mail : calibratech-cal@yahoo.com, calibratech-cal@hotmail.com

Certificate of Calibration

Certificate No. : 64-420032-1

Page : 2 of 2

Result of Calibration :

UUC Condition As-Received : Good

Function : Electrical measurement

pH meter

Performing standard curve by Multiproduct Calibrator at pH (4,7,10)

Adjustment Curve at nominal pH	Applied Voltage (mV)	Nominal Value (pH)	UUC Reading		Correction (mV)	Uncertainty (± mV)
			(pH)	(mV)		
4, 7, 10	177.4800	4	4.00	177	0	0.58
	0.0000	7	7.00	0	0	0.58
	-177.4800	10	10.00	-178	1	0.58

Function : pH meter with electrode

Performing a three - buffer standard curve using buffer nominal pH (4,7,10)

Adjustment Curve at nominal pH	Standard Buffer (pH)	UUC Reading (pH)	Correction (pH)	Uncertainty (± pH)
4, 7, 10	4.004	4.00	0.00	0.011
	6.985	7.00	-0.01	0.020
	9.963	10.00	-0.04	0.053

Remark

UUC : Unit Under Calibration

This result of calibration was found accurate as shown on date and place of calibration only.

This reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$,
providing a level of confidence of approximately 95%

- 000 -





National Food Institute, Ministry of Industry, Thailand

2008 Soi Arun Ammarin 36, Arun Ammarin Rd., Bangyeekhan, Bangphlad, Bangkok 10700 Thailand
Tel : +66 (0) 2422 8688 Fax : +66 (0) 2422 8558 Website : www.nfi.or.th E-mail : cal@nfi.or.th



Calibration Certificate

Certificate No.: 2004097-002-01
Client name: SECOT CO., LTD.
Address: 239 Rimklongprapa Road, Bangsue,
Bangsue, Bangkok 10800

Page 1 of 3

Equipment: CHAMBER (Incubator)

Manufacturer: MEMMERT

Model: ICP 400

Serial No.: K406.0004

ID No.: N/A

Order No.: 2004097

Operation No.: 2004097-002

Date of Receipt: 3 August 2020

Date of Calibration: 3 August 2020

Calibrated by Mr.Pheraphat Tuanjit
Scientist

Approved by


(Mr.Manas Somsak)

Senior Analyst, Division of Calibration Laboratory

Date of Issue: 5 August 2020

Responsible for the Technical Management Team

The uncertainties are for a confidence probability of approximately 95 %.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.



National Food Institute, Ministry of Industry, Thailand

2008 Soi Arun Ammarin 36, Arun Ammarin Rd., Bangyeekhan, Bangphlad, Bangkok 10700 Thailand
Tel : +66 (0) 2422 8688 Fax : +66 (0) 2422 8558 Website : www.nfi.or.th E-mail : cal@nfi.or.th



Calibration Report

Certificate No.: 2004097-002-01
Equipment: CHAMBER (Incubator)
Model: ICP 400 Serial No.: K406.0004
Resolution: 0.1 °C ID No.: N/A
Manufacturer: MEMMERT

Date of Calibration: 3 August 2020

Page 2 of 3

Location: Walkway, SECOT CO., LTD.
Environment Condition: Ambient Temperature (31 ± 1) °C
Relative Humidity (55 ± 5) %
Line Voltage (229 ± 1) Volt

Condition of this results of Calibration:

- This instrument was calibrated by insert 9 standard thermometer into its chamber and calibration according to W-TE-014 Based on TLAS G-20-1/02-08 (E); Guidelines for Calibration and Checks of Temperature Controlled Enclosures.
 - The temperature scale used was based on ITS - 90.
 - All data show below were final values and the initial data may be obtained upon request.

2. Reference Standard Instrument :

Instrument	Model	Serial No./ID No.	Certificate No.	Due Date	Through
Digital Thermometer with sensor	34972A	MY49020132	TE 630168-01	27 January 2021	NATIONAL FOOD INSTITUTE
	RTD	CH#101-109/ RTD#101-109			

- This certificate is traceable to International System of Units (SI Units).
- This certificate was certified only for the instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.
- Condition of Calibrated item : Good

UUC Description :

Time of Record 1 Hour 9 Minute At 20.0 °C
Fresh air Damper ☐ Open Position ☐
☒ Close
☐ Not Available

7. Result of Calibration : ☒ Without adjustment ☐ After adjustment

Am -

Calibration Report

Certificate No.: 2004097-002-01
Equipment: CHAMBER (Incubator)
Model: ICP 400 Serial No.: K406.0004
Resolution: 0.1 °C ID No.: N/A
Manufacturer: MEMMERT

Date of Calibration: 3 August 2020

Page 3 of 3

Calibration point: 20.0 °C

Calibration result:

Calibration Condition	Temperature (°C)	Relative Humidity (%)	Line Voltage (Volt)
MIN	31.1	50	228.8
MAX	31.6	60	230.0

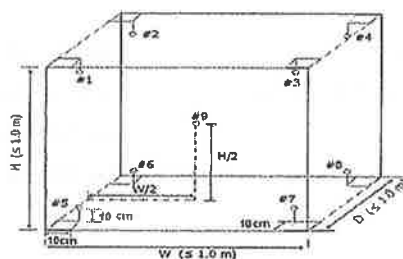


Table 1 : Reporting of Temperature

Calibration point (°C)	Measured Temperature (°C) @ Sensor No. (Sensor No.9 is REF)									Uncertainty ± (°C)
	# 1	# 2	# 3	# 4	# 5	# 6	# 7	# 8	# 9	
20.0	19.94	20.08	20.14	20.09	20.13	20.11	20.21	20.01	20.08	0.27

Table 2 : Reporting of Characterization Result

UUC* Setting (°C)	UUC* reading (°C)			Stability ± (°C)	Uniformity (°C)	Overall Variation (°C)
	MIN	MAX	Average			
20.0	19.9	20.0	20.0	0.06	0.16	0.38

Note The quoted uncertainty include " Stability " and " Loading effect (20% of Temp Uniformity) "

UUC* = Unit Under Calibration

Stability = One-half of the greatest maximum difference of measured temperatures at any one sensors, for at least half an hour after reaching steady state.

Uniformity = The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time.

Overall Variation = The difference of the maximum and minimum measured temperatures throughout observation time.

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor $k=2$, providing a level of confidence of approximately 95 %.

----- End -----

Signature



National Food Institute, Ministry of Industry, Thailand

2008 Sol Arun Ammarin 36, Arun Ammarin Rd., Bangyeekhan, Bangphlad, Bangkok 10700 Thailand
Tel : +66 (0) 2422 8688 Fax : +66 (0) 2422 8558 Website : www.nfi.or.th E-mail : cal@nfi.or.th



Calibration Certificate


Certificate No.: 2004097-003-01
Client name: SECOT CO., LTD.
Address: 239 Rimklongprapa Road, Bangsue,
Bangsue, Bangkok 10800

Page 1 of 3

Equipment: Water Bath
Manufacturer: MEMMERT
Model: WB 29
Serial No.: 1698.0051
ID No.: N/A
Order No.: 2004097
Operation No.: 2004097-003
Date of Receipt: 3 August 2020
Date of Calibration: 3 August 2020

Calibrated by Mr.Pheraphat Tuanjit
Scientist

Approved by


(Mr.Manas Somsak)

Senior Analyst, Division of Calibration Laboratory
Responsible for the Technical Management Team

Date of Issue: 5 August 2020

The uncertainties are for a confidence probability of approximately 95 %.

This Certificate is issued in accordance with the conditions of accreditation granted by the National Accreditation System of Thailand which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.



National Food Institute, Ministry of Industry, Thailand

2008 Soi Arun Ammarin 36, Arun Ammarin Rd., Bangyeekhan, Bangphlad, Bangkok 10700 Thailand
Tel : +66 (0) 2422 8688 Fax : +66 (0) 2422 8558 Website : www.nfi.or.th E-mail : cal@nfi.or.th



Calibration Report

Certificate No.: 2004097-003-01
Equipment: Water Bath
Model: WB 29 Serial No.: 1698.0051
Resolution: 0.1 °C ID No.: N/A
Manufacturer: MEMMERT
Date of Calibration: 3 August 2020

Page 2 of 3

Location: Walkway, SECOT CO., LTD.
Environment Condition: Ambient Temperature (32 ± 1) °C
Relative Humidity (45 ± 5) %
Line Voltage (227 ± 2) Volt

Condition of this results of Calibration:

- This instrument was calibrated by insert 5 standard thermometer into its liquid bath and calibration according to W-TE-011 based on ASTM E715-80 (Reapproved 2006): Standard Specification for Gravity-Convection and Forced-Circulation Water Baths.
 - The temperature scale used is ITS - 90.
 - All data show below were final values and the initial data may be obtained upon request.

2. Reference Standard Instrument :

Instrument	Model	Serial No./ID No.	Certificate No.	Due Date	Through
Digital Thermometer with sensor	34972A	MY49020132	TE 630168-01	27 January 2021	NATIONAL FOOD INSTITUTE
	RTD	RTD#201-205 / CH#201-205			

- This certificate is traceable to International System of Units (SI Units).
- This certificate was certified only for the instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.
- Condition of Calibrated item : Good

UUC Description:

Time of Record 1 Hour 9 Minute At 95.0 °C

7. Result of Calibration :
- | | |
|-------------------------------------|--------------------|
| <input checked="" type="checkbox"/> | Without adjustment |
| <input type="checkbox"/> | After adjustment |

Signature

Calibration Report

Certificate No.: 2004097-003-01

Equipment: Water Bath

Model: WB 29

Serial No.: 1698.0051

Resolution: 0.1 °C

ID No.: N/A

Manufacturer: MEMMERT

Date of Calibration: 3 August 2020

Page 3 of 3

Calibration point: 95.0 °C

Calibration result:

calibration Condition	Temperature (°C)	Relative Humidity (%)	Line Voltage (Volt)
Min	31.5	40	225.0
Max	32.2	50	228.0

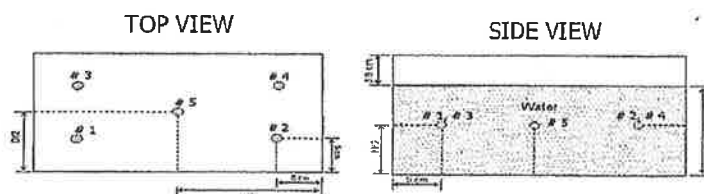


Table 1 : Reporting of Temperature

Sensor Installation Location

Calibration Point (°C)	Measured Temperature (°C) @ Sensor No. (Sensor No.5 is REF)					Uncertainty ± (°C)
	# 1	# 2	# 3	# 4	# 5	
95.0	95.19	94.93	95.01	94.99	95.12	0.34

Table 2 : Reporting of Characterization Result

UUC* Setting (°C)	UUC* reading (°C)			Stability ± (°C)	Uniformity (°C)	Overall Variation (°C)
	MIN	MAX	Average			
95.0	94.9	95.2	95.0	0.20	0.33	0.55

Note The quoted uncertainty include " Stability " and " Loading effect (20% of Temp Uniformity)"

UUC* = Unit Under Calibration

Stability = One-half of the greatest maximum difference of measured temperatures at any one sensors,
for at least half an hour after reaching steady state.

Uniformity = The maximum difference of measured temperatures at any sensors and the measured
temperature at the reference location which are observed at the same time.

Overall Variation = The difference of the maximum and minimum measured temperatures throughout observation time.

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k= 2,
providing a level of confidence of approximately 95 %.

----- End -----

**SOUND LEVEL METER CALIBRATION**

Calibration Location: SECOT

Calibration Date: Feb 18, 21

SOUND LEVEL CALIBRATOR

Brand	Model	Serial No.	Calibrated (dB)	Frequency (Hz)		
RION	NC-74	34283648	94.00	1000		

No.	Brand	Model	Serial No.	Microphone Serial No.	SLM Reading (dB)	dB Adjust
71	RION	NL-21	00487728	119000	94.2	-0.2
92	RION	NL-21	00198274	123477	94.2	-0.2

Calibrated by :

Approved by :

Preeda S.



**ELECTRICAL AND ELECTRONICS INSTITUTE
FOUNDATION FOR INDUSTRIAL DEVELOPMENT**

975 Moo 4, Bangpoo Industrial Estate, Soi 8, Sukhumvit Road km 37,

Phraek Sa, Mueang Samut Prakan, Samut Prakan 10280

Tel: +66 2709 4860-8 Fax: +66 2324 0917-8



Certificate No.: 0037SV21
Operation No.: CP2021010002

Certificate of Calibration

Equipment: Sound Calibrator
Manufacturer: RION
Model/Type: NC-74
Serial No.: 34283648
ID No.: -
Customer: SECOT Co.,Ltd.
Address: 239 Rimklongprapa Rd., Bangsue,
Bangkok 10800 Thailand
Received Date: 4 January 2021
Calibrated Date: 6 January 2021
Issued Date: 8 January 2021
Calibrated by: Ms. Juntaporn Kunhakom

Approved by: _____

(Mr. Sittichai Swaksuriya-wong)

Group Manager
สถาบันไฟฟ้าและอิเล็กทรอนิกส์
ELECTRICAL AND ELECTRONICS INSTITUTE

The reported uncertainty of measurement was based on standard uncertainty multiplied by a coverage factor $k = 2.00$, providing a level of confidence of approximately 95%. This certificate may not be reproduced other than in full except with the prior written approval of the Electrical and Electronics Institute, Foundation for Industrial Development.

Certificate No.: 0037SV21

Calibration Report

Equipment: Sound Calibrator
Manufacturer: RION
Model/Type: NC-74
Serial No.: 34283648
ID No.: -
Ambient Temperature: (23 ± 2) °C
Relative Humidity: (50 ± 15) %
Pressure: (101.3 ± 1.5) kPa
Method of Calibration :-
IEC 60942:2017

Condition of this result of calibration

1. Reference standards instrument :-

<u>Instrument</u>	<u>Model</u>	<u>Serial No.</u>	<u>Cert. No.</u>	<u>Due Date</u>
1) Standard microphone	4180	2661000	AA-1013-20	12 May 2021
2) Waveform Generator	33511B	MY52302264	0100RF20	17 June 2021
3) Audio Analyzing DMM	2015-P	000136E	E1U203927	16 November 2021
4) Pressure humidity and Temperature Transmitter	PTU301	F0640002	CL1-P200051 0305TE20	31 May 2021 29 June 2021

2. This result of calibration was found accurate as shown on date and place of calibration only.

3. This certification is traceable to the international system of unit maintained at :-

Reference standards instrument for Acoustic function

- National Institute of Metrology (Thailand)

Reference standards instrument for Electrical function

- Electrical and Electronics Institute; ONSC Accredited Calibration No.0119

Result of Calibration:-

1. Function : Sound pressure level

Norminal Frequency (Hz)	Specified Sound Pressure level (dB)	Measured value (dB)	Deviated value ^[1] (dB)	Acceptance limit ^[3] (dB)
1000	94	94.24	0.24	±0.25

2. Function : Frequency

Norminal Sound Pressure level (dB)	Specified Frequency (Hz)	Measured value (Hz)	Deviated value ^[2] (%)	Acceptance limit ^[3] (%)
94	1000	1003.2	0.3	±0.7

Certificate No.: 0037SV21

Calibration Report

3. Function : Total distortion + noise

Normal	Normal	Measured value ^[4]	Acceptance limit ^[5]
Sound Pressure level (dB)	Frequency (Hz)	(%)	(%)
94	1000	1.3	2.5

Uncertainty of measurement

Function	Uncertainty	Maximum-permitted uncertainty of measurement
Sound pressure level	0.10 dB	0.15 dB
Frequency	0.10 %	0.20 %
Total distortion + noise	0.40 %	0.50 %

- Note:
- [1] The deviated value is the absolute value of the difference between the measured value and the corresponding specified sound pressure level.
 - [2] The deviated value is the absolute value of the difference in percent between the measured value and the corresponding specified frequency.
 - [3] The acceptance limit is for the deviated value.
 - [4] The measured value is the total distortion + noise, measured over the frequency range from 20 Hz to 20 kHz.
 - [5] The acceptance limit is for the Measured value.

Remarks: 1. Using the 1/2-inch microphone adaptor NC-74-002.
2. Acceptance limit was IEC 60942:2017 Class 1.

-- End of Report --



SOUND LEVEL METER CALIBRATION

Calibration Location: **SECOT**

Calibration Date: **Feb 18, 21**

SOUND LEVEL CALIBRATOR

Brand	Model	Serial No.	Calibrated (dB)	Frequency (Hz)		
CASELLA	CEL120/2	2839225	114.0	1000		
No.	Brand	Model	Serial No.	Microphone Serial No.	SLM Reading (dB)	dB Adjust
22	CASELLA	CEL-246	3173338	3173338	114.0	0.0
23	CASELLA	CEL-246	3173339	3173339	113.9	0.1

Calibrated by :

Approved by :

Certificate No.: 0038SV21
Operation No.: CP2021010003

Certificate of Calibration

Equipment: Sound Calibrator
Manufacturer: CASELLA
Model/Type: CEL-120/2
Serial No.: 2839225
ID No.: -
Customer: SECOT Co.,Ltd.
Address: 239 Rimklongprapa Rd., Bangsue,
Bangkok 10800 Thailand
Received Date: 4 January 2021
Calibrated Date: 6 January 2021
Issued Date: 8 January 2021
Calibrated by: Ms. Juntaporn Kunhakom

Approved by:

(Mr. Sittichai Swaksuriyawong)

Group Manager
สถาบันไฟฟ้าและอิเล็กทรอนิกส์
ELECTRICAL AND ELECTRONICS INSTITUTE

The reported uncertainty of measurement was based on standard uncertainty multiplied by a coverage factor $k = 2.00$, providing a level of confidence of approximately 95%. This certificate may not be reproduced other than in full except with the prior written approval of the Electrical and Electronics Institute, Foundation for Industrial Development.

Certificate No.: 0038SV21

Calibration Report

Equipment: Sound Calibrator
Manufacturer: CASELLA
Model/Type: CEL-120/2
Serial No.: 2839225
ID No.: -
Ambient Temperature: (23 ± 2) °C
Relative Humidity: (50 ± 15) %
Pressure: (101.3 ± 1.5) kPa

Method of Calibration :-

IEC 60942:2017

Condition of this result of calibration

1. Reference standards instrument :-

Instrument	Model	Serial No.	Cert. No.	Due Date
1) Standard microphone	4180	2661000	AA-1013-20	12 May 2021
2) Waveform Generator	33511B	MY52302264	0100RF20	17 June 2021
3) Audio Analyzing DMM	2015-P	000136E	E1U203927	16 November 2021
4) Pressure humidity and Temperature Transmitter	PTU301	F0640002	CL1-P200051 0305TE20	31 May 2021 29 June 2021

2. This result of calibration was found accurate as shown on date and place of calibration only.

3. This certification is traceable to the international system of unit maintained at :-

Reference standards instrument for Acoustic function

- National Institute of Metrology (Thailand)

Reference standards instrument for Electrical function

- Electrical and Electronics Institute; ONSC Accredited Calibration No.0119

Result of Calibration:-

1. Function : Sound pressure level

Normalinal Frequency (Hz)	Specified Sound Pressure level (dB)	Measured value (dB)	Deviated value ^[1] (dB)	Acceptance limit ^[3] (dB)
1000	114	114.17	0.17	±0.40

2. Function : Frequency

Normalinal Sound Pressure level (dB)	Specified Frequency (Hz)	Measured value (Hz)	Deviated value ^[2] (%)	Acceptance limit ^[3] (%)
114	1000	1000.0	0.0	±1.7

Certificate No.: 0038SV21

Calibration Report

3. Function : Total distortion + noise

Normal Sound Pressure level (dB)	Normal Frequency (Hz)	Measured value ^[4] (%)	Acceptance limit ^[5] (%)
114	1000	0.3	3.0

Uncertainty of measurement

Function	Uncertainty	Maximum-permitted uncertainty of measurement
Sound pressure level	0.10 dB	0.35 dB
Frequency	0.10 %	0.20 %
Total distortion + noise	0.40 %	1.00 %

- Note:
- [1] The deviated value is the absolute value of the difference between the measured value and the corresponding specified sound pressure level.
 - [2] The deviated value is the absolute value of the difference in percent between the measured value and the corresponding specified frequency.
 - [3] The acceptance limit is for the deviated value.
 - [4] The measured value is the total distortion + noise, measured over the frequency range from 20 Hz to 20 kHz.
 - [5] The acceptance limit is for the Measured value.

Remarks: 1. Acceptance limit was IEC 60942:2017 Class 2.

-- End of Report --



Heat Stress Meter Calibration

Date: Mar 26,21

Temperature (°C) 20

Barometric Pressure: Pb (mmHg) 760

REFERENCE STANDARD INSTRUMENT

UNIT UNDER TEST

Equipment : DIGITAL THERMOCOUPLE CALIBRATOR

Equipment : Heat Stress Meter

Model No. 714

Model No. QUESTEMP °46

Serial No. 7590122

Serial No. TSR010003

Manufacturer FLUKE

Manufacturer 3M/QUEST Technology

Calibration Date 15 July 2020

Temperature Reading

Reference Setting (°C)	Tg (°C)	T (°C)	Tn (°C)
30.0	30.2	30.2	30.1
40.0	40.1	40.3	40.1
50.0	50.3	50.4	50.2
60.0	60.2	60.2	60.3

- Note : 1) Tg = Globe thermometer temperature
 2) Tn = Wet bulb with natural ventilation temperature
 3) T = Ambient temperature

Calibrated by :

Approved by :



Heat Stress Meter Calibration

Date: Mar 25,21

Temperature (°C) 20

Barometric Pressure: Pb (mmHg) 760

REFERENCE STANDARD INSTRUMENT

UNIT UNDER TEST

Equipment : DIGITAL THERMOCOUPLE CALIBRATOR

Equipment : Heat Stress Meter

Model No. 714

Model No. QUESTEMP °46

Serial No. 7590122

Serial No. TSM050001

Manufacturer FLUKE

Manufacturer 3M/QUEST Technology

Calibration Date 15 July 2020

Temperature Reading

Reference Setting (°C)	T _g (°C)	T (°C)	T _n (°C)
30.0	30.3	30.1	30.2
40.0	39.9	40.0	39.9
50.0	49.7	49.9	49.9
60.0	59.8	59.8	59.8

- Note : 1) T_g = Globe thermometer temperature
 2) T_n = Wet bulb with natural ventilation temperature
 3) T = Ambient temperature

Calibrated by : Suraphong P.

Approved by : Wittaya K.



Heat Stress Meter Calibration

Date: Mar 25,21

Temperature (°C) 20

Barometric Pressure: Pb (mmHg) 760

REFERENCE STANDARD INSTRUMENT

UNIT UNDER TEST

Equipment : DIGITAL THERMOCOUPLE CALIBRATOR

Equipment : Heat Stress Meter

Model No. 714

Model No. QUESTEMP °46

Serial No. 7590122

Serial No. TSM050004

Manufacturer FLUKE

Manufacturer 3M/QUEST Technology

Calibration Date 15 July 2020

Temperature Reading

Reference Setting (°C)	Tg (°C)	T (°C)	Tn (°C)
30.0	30.1	29.9	30.0
40.0	40.2	40.2	40.1
50.0	50.2	50.3	50.2
60.0	60.2	60.2	60.3

- Note : 1) Tg = Globe thermometer temperature
 2) Tn = Wet bulb with natural ventilation temperature
 3) T = Ambient temperature

Calibrated by :

Approved by :



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG, BANGKOK 10250
TEL. 0-2717-3000-24 FAX. 0-2719-9484



Certificate of Calibration

Certificate No. : 20T1520

Page : 1 of 2

Equipment : Digital Thermometer With Sensor

Manufacturer: Fluke

Model : 714

Serial No.: 7590122

ID No.: 200209-1

Condition As-Received: Used Item

Received Date: 24 June 2020

Calibration Date: 13 July 2020
to 15 July 2020

Reference: 2006-0946DN

Ambient Temperature: (25 ± 3) °C

Relative Humidity: (50 ± 20) %

This certificate may not be reproduced other than in full,
except with the prior written approval of the head of
Corporate Services 3: Equipment Calibration and Testing Services.

Submitted by: Secot Co.,Ltd.

239 Rimklongprapa Road., Bangsue, Bangkok 10800

Procedure used: Calibration were conducted using in-house calibration procedure CP-T01 according to comparison with Industrial Platinum Resistance Thermometer (IPRT) into liquid bath temperature controller.
The temperature scale used was based on ITS-90.

Condition of this result of calibration

1.Reference standards instruments :

<u>Instrument</u>	<u>Model</u>	<u>Serial No.</u>	<u>Certificate No.</u>	<u>Due Date</u>
1) Digital Thermometer	1529	A7A609	19I1235	24 Sep 2020
2) Industrial Platinum Resistance Thermometer	5627-12	571975	19I1235	24 Sep 2020
3) Industrial Platinum Resistance Thermometer	5627-12	571970	19I1235	24 Sep 2020
4) Industrial Platinum Resistance Thermometer	5627-12	555541	19I1235	24 Sep 2020
5) Digital Thermometer	1529	A66176	19I1397	01 Nov 2020
6) Industrial Platinum Resistance Thermometer	5627	739437	19I1397	01 Nov 2020

2.The certificate is valid only to the item calibrated on date and place of calibration.

3.This Certification is traceable to the International System of Unit maintained at:-

-National Institute of Metrology Thailand (NIMT)

Calibrated by : Yossapon Poljorn

Issue Date : 20 July 2020

Approved Signatory :

[] Phalinee Prabpaipal

[] Chatchawan Khunpiluek

[✓] Wanlop Larpkurn

B 0234488



Cert. No.: 20T1520

Page.: 2 of 2

Result of Calibration:-

Without Adjustment

Function:

Temperature measurement

This equipment was connected with Thermocouple Type K

ID No. 200209-1

Immersion	Standard	UUC*		Uncertainty
Depth	Temperature	Reading	Error	of Measurement
(mm.)	(°C)	(°C)	(°C)	(±°C)
150	-0.0007	0.7	0.7007	0.38
150	20.0056	20.1	0.0944	0.39
150	80.0019	79.0	-1.0019	0.43
150	104.0042	103.1	-0.9042	0.50
150	110.0051	108.9	-1.1051	0.52
150	150.0041	149.8	-0.2041	0.63
150	180.0045	180.2	0.1955	0.72

UUC* : Unit Under Calibration

The reported uncertainty of measurement was based on standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%.

-o0o-

ภาคผนวก ก

หนังสือรับรองผู้รับรองรายงานการตรวจวัดและวิเคราะห์
สถานะการทำงาน ตามกฎกระทรวง กำหนดมาตรฐานในการบริหาร
และการจัดการด้านความปลอดภัย อาชีวอนามัย และสภาพแวดล้อม
ในการทำงาน เกี่ยวกับความร้อน แสงสว่าง และเสียง พ.ศ.2549



เลขที่ ๑๐๗/๒๕๕๗




กรมสวัสดิการและคุ้มครองแรงงาน
ถนนมิตรไมตรี เขตดินแดง
กรุงเทพมหานคร ๑๐๔๐๐

หนังสือฉบับนี้ให้ไว้เพื่อรับรองว่านางสาวสุนันทา ศิริพัฒน์ นาม น. เลขประจำตัวประชาชน ๓ ๑๐๐๘ ๐๐๖๓๗ ๒๗ ๙ เลขทะเบียน รสส. ๐๐๔-๕๗/๐๔๒๙ เป็นผู้ที่สามารถรับรองรายงาน การตรวจวัดและวิเคราะห์สภาวะการทำงานตามกฎกระทรวงกำหนดมาตรฐานในการบริหารและการจัดการ ด้านความปลอดภัย อาชีวอนามัย และสภาพแวดล้อมในการทำงานเกี่ยวกับความร้อน แสงสว่าง และเสียง พ.ศ. ๒๕๕๙ ลงวันที่ ๑๖ กุมภาพันธ์ พ.ศ. ๒๕๕๙ หมวด ๕ ข้อ ๑๖ ซึ่งมีผลใช้บังคับตามพระราชบัญญัติ ความปลอดภัย อาชีวอนามัย และสภาพแวดล้อมในการทำงาน พ.ศ. ๒๕๕๔ มาตรา ๗๔ ตั้งแต่วันที่ ๓๐ ตุลาคม พ.ศ. ๒๕๕๗ ถึงวันที่ ๒๙ ตุลาคม พ.ศ. ๒๕๕๘

ให้ไว้ ณ วันที่ ๓๐ ตุลาคม พ.ศ. ๒๕๕๗



(นายวรณัฐ ปิตุวรรณ)
รองอธิบดี ปฏิบัติราชการแทน
อธิบดีกรมสวัสดิการและคุ้มครองแรงงาน


(นางสาวสุนันทา ศิริพัฒน์ นาม น.)

เลขที่ ๐๙๖/๒๕๕๘

กรมสวัสดิการและคุ้มครองแรงงาน
ถนนมิตรไมตรี เขตดินแดง
กรุงเทพมหานคร ๑๐๔๐๐

หนังสือฉบับนี้ให้ไว้เพื่อรับรองว่านางสาวสุนันทา ศิริพัฒน์ นาม น. เลขประจำตัวประชาชน ๓ ๑๐๐๘ ๐๐๖๓๗ ๒๗ ๙ เลขทะเบียน รสส. ๐๐๔-๕๗/๐๔๒๙ เป็นผู้ที่สามารถรับรองรายงาน การตรวจวัดและวิเคราะห์สภาวะการทำงานตามกฎกระทรวงกำหนดมาตรฐานในการบริหารและการจัดการ ด้านความปลอดภัย อาชีวอนามัย และสภาพแวดล้อมในการทำงานเกี่ยวกับความร้อน แสงสว่าง และเสียง พ.ศ. ๒๕๕๙ ลงวันที่ ๑๖ กุมภาพันธ์ พ.ศ. ๒๕๕๙ หมวด ๕ ข้อ ๑๖ ซึ่งมีผลใช้บังคับตามพระราชบัญญัติ ความปลอดภัย อาชีวอนามัย และสภาพแวดล้อมในการทำงาน พ.ศ. ๒๕๕๔ มาตรา ๗๔ ตั้งแต่วันที่ ๓๐ ตุลาคม พ.ศ. ๒๕๕๗ ถึงวันที่ ๒๙ ตุลาคม พ.ศ. ๒๕๕๘

ให้ไว้ ณ วันที่ ๒๙ กันยายน พ.ศ. ๒๕๕๘



(นายวรณัฐ ปิตุวรรณ)
รองอธิบดี ปฏิบัติราชการแทน
อธิบดีกรมสวัสดิการและคุ้มครองแรงงาน


(นางสาวสุนันทา ศิริพัฒน์ นาม น.)



แบบ รสส. ๒

เลขที่ ๑๒๖/๒๕๕๙

กรมสวัสดิการและคุ้มครองแรงงาน
ถนนมิตรไมตรี เขตดินแดง
กรุงเทพมหานคร ๑๐๔๐๐

หนังสือฉบับนี้ให้ไว้เพื่อรับรองว่านางสาวสุนันทา ศิริภูมินานนท์ เลขประจำตัวประชาชน ๓ ๑๐๐๘ ๐๐๖๓๗ ๒๗ ๙ เลขทะเบียน รสส. ๐๐๔-๕๙/๐๓๒๕ เป็นผู้ที่สามารถรับรองรายงานการตรวจวัดและวิเคราะห์สภาวะการทำงานตามกฎหมายกระทรวงกำหนดมาตรฐานในการบริหารและการจัดการด้านความปลอดภัย อาชีวอนามัย และสภาพแวดล้อมในการทำงานเกี่ยวกับความร้อน แสงสว่าง และเสียง พ.ศ. ๒๕๕๙ ลงวันที่ ๑๖ กุมภาพันธ์ พ.ศ. ๒๕๕๙ หมวด ๕ ข้อ ๑๖ ซึ่งมีผลใช้บังคับตามพระราชบัญญัติความปลอดภัย อาชีวอนามัย และสภาพแวดล้อมในการทำงาน พ.ศ. ๒๕๕๙ มาตรา ๗๔ ตั้งแต่วันที่ ๓๐ ตุลาคม พ.ศ. ๒๕๕๙ ถึงวันที่ ๒๙ ตุลาคม พ.ศ. ๒๕๖๐

ให้ไว้ ณ วันที่ ๑๐ ตุลาคม พ.ศ. ๒๕๕๙

(นายอำนาจ ภูระหงษ์)
ผู้อำนวยการกองความปลอดภัยแรงงาน



(นางสาวสุนันทา ศิริภูมินานนท์)

ภาคผนวก ข

หนังสืออนุญาตขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
จากกรมโรงงานอุตสาหกรรม



ที่อก ๐๓๑๐(๑)/ ๑๑ ๘ ๐ ๕

กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ เขตราชเทวี
กรุงเทพมหานคร ๑๐๔๐๐

๒๑ ตุลาคม ๒๕๖๓

เรื่อง ค่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

เรียน กรรมการผู้จัดการ บริษัท ซิคอท จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๒๔ มีนาคม ๒๕๖๓

- สิ่งที่ส่งมาด้วย ๑. รายชื่อผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๑ แผ่น
๒. รายชื่อเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๑ แผ่น
๓. ขอบข่ายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๓๒ แผ่น

ตามหนังสืออ้างถึง บริษัท ซิคอท จำกัด ค่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ๖-๒๓๙ สถานที่ตั้งเลขที่ ๒๓๙ ถนนริมคลองประปา แขวงบางซื่อ เขตบางซื่อ กรุงเทพมหานคร ต่อกรมโรงงานอุตสาหกรรม นั้น


กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้บริษัท ซิคอท จำกัด ค่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน โดยมีองค์ประกอบดังนี้

- ก. ผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๑๐ ราย ตามสิ่งที่ส่งมาด้วย ๑
ข. เจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๓๑ ราย ตามสิ่งที่ส่งมาด้วย ๒
ค. ขอบข่ายสารมลพิษที่ได้รับขึ้นทะเบียนให้วิเคราะห์ในน้ำเสีย จำนวน ๔๖ รายการ น้ำได้ดิน จำนวน ๑๒๓ รายการ อากาศเสีย จำนวน ๒๗ รายการ สิ่งปฏิกูลหรือวัสดุที่ไม่ใช้แล้ว จำนวน ๓๔ รายการ และดิน จำนวน ๑๒๒ รายการ รวมทั้งสิ้นจำนวน ๓๕๒ รายการ ตามสิ่งที่ส่งมาด้วย ๓

หนังสือฉบับนี้จะหมดอายุในวันที่ ๒ พฤษภาคม ๒๕๖๖ หากประสงค์จะต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ให้ยื่นคำขอต่ออายุพร้อมเอกสารประกอบคำขอต่อกรมโรงงานอุตสาหกรรมภายใน ๓๐ วัน ก่อนวันสิ้นอายุของหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ซึ่งคำขอต่ออายุดังกล่าวขอรับได้ที่กรมโรงงานอุตสาหกรรม

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

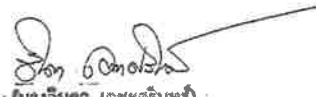
(นางจันทา เดชะธื่นท)
ผู้อำนวยการกองวิจัยและเฝ้าระวังมลพิษโรงงาน
ปฏิบัติการตามพันธกิจกรมโรงงานอุตสาหกรรม

กองวิจัยและเฝ้าระวังมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๒๐๒ ๔๐๐๒ ๐ ๒๒๐๒ ๔๑๔๖

โทรสาร ๐ ๒๓๕๕ ๓๒๐๘ ๐ ๒๓๕๕ ๓๔๑๕

เอกสารแนบท้ายหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท ซิคอท จำกัด

เลขทะเบียน ๖-๒๓๙

ที่อก ๐๓๑๐(๑)/ ๑๑ ๘ ๐ ๕

ลงวันที่ ๒๑ ตุลาคม ๒๕๖๓

ก. ผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๑๐ ราย

- | | |
|-----------------------------|----------------------------|
| ๑) นางสาวสุดี เกรียงไกรอุดม | ทะเบียนเลขที่ ๖-๒๓๙-ค-๕๘๒๐ |
| ๒) นางอารยา ทิพรักษ์ | ทะเบียนเลขที่ ๖-๒๓๙-ค-๕๘๖๓ |
| ๓) นายชรรชัย เกรียงไกรอุดม | ทะเบียนเลขที่ ๖-๒๓๙-ค-๕๘๗๕ |
| ๔) นางสาวเชมชุตา อินทร์ศรี | ทะเบียนเลขที่ ๖-๒๓๙-ค-๕๘๗๖ |
| ๕) นางสาวปรีดา สมใจ | ทะเบียนเลขที่ ๖-๒๓๙-ค-๕๘๗๘ |
| ๖) นางสาวอรัญญา มาตา | ทะเบียนเลขที่ ๖-๒๓๙-ค-๕๘๗๙ |
| ๗) นางสาวลดาวัลย์ วงศ์เจริญ | ทะเบียนเลขที่ ๖-๒๓๙-ค-๕๘๘๐ |
| ๘) นางสาวมณีนรณ เกตวันดี | ทะเบียนเลขที่ ๖-๒๓๙-ค-๕๘๘๒ |
| ๙) นางสาวนริสา ภูวสรเพ็ญ | ทะเบียนเลขที่ ๖-๒๓๙-ค-๖๔๑๙ |
| ๑๐) นางสาวศิริวรรณ ฉิมสง่า | ทะเบียนเลขที่ ๖-๒๓๙-ค-๖๔๒๐ |

เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท ซีคอฟ จำกัด

เลขทะเบียน ว-๒๓๙

ที่ อก ๐๓๓๐(๑)/ ๕ ๑ ๘ ๐ ๔

ลงวันที่ ๒๑ ตุลาคม ๒๕๖๓

ข. เจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๓๑ ราย

๑) นางสาวสุวิชาติ ชัยธรรม	ทะเบียนเลขที่	ว-๒๓๙-จ-๕๕๒๕
๒) นางสาวสุธาทิพย์ เทียนเตี้ย	ทะเบียนเลขที่	ว-๒๓๙-จ-๕๕๒๙
๓) นางสาวสุนันดา ศิริวัฒนานนท์	ทะเบียนเลขที่	ว-๒๓๙-จ-๕๕๔๓
๔) นายบวร ศิชัยยะ	ทะเบียนเลขที่	ว-๒๓๙-จ-๕๕๔๖
๕) นางสาวเกศรินทร์ วรเดชาวิทยา	ทะเบียนเลขที่	ว-๒๓๙-จ-๕๕๔๑
๖) นายอนันต์พนัน พินวันนา	ทะเบียนเลขที่	ว-๒๓๙-จ-๖๐๐๑
๗) นายชิตพล สมประสงค์	ทะเบียนเลขที่	ว-๒๓๙-จ-๖๐๐๒
๘) นางสาวศศิธร พรหมประเสริฐ	ทะเบียนเลขที่	ว-๒๓๙-จ-๖๐๐๓
๙) นายศิวะนนท์ กุลวงษ์	ทะเบียนเลขที่	ว-๒๓๙-จ-๖๐๐๕
๑๐) นางสาวโชติมาส ไทยเจริญ	ทะเบียนเลขที่	ว-๒๓๙-จ-๖๐๐๖
๑๑) นางสาวปิยขวัญ สุระโคตร	ทะเบียนเลขที่	ว-๒๓๙-จ-๖๔๒๑
๑๒) นางสาวณัฐศิริ เลิศธีรพิพัฒน์	ทะเบียนเลขที่	ว-๒๓๙-จ-๖๔๒๓
๑๓) นางสาวเกษวรินทร์ ศิลศึก	ทะเบียนเลขที่	ว-๒๓๙-จ-๖๔๒๔
๑๔) นางสาวอลิษา คนิวรานนท์	ทะเบียนเลขที่	ว-๒๓๙-จ-๗๒๓๑
๑๕) นางสาวจิรณันท์ จิตตะศรี ปิยะธนากร	ทะเบียนเลขที่	ว-๒๓๙-จ-๗๒๓๒
๑๖) นางสาวสิริวรรณ แก้วชิงดวง	ทะเบียนเลขที่	ว-๒๓๙-จ-๗๒๓๓
๑๗) นางสาวปัทมวรรณ สุวรรณโรจน์	ทะเบียนเลขที่	ว-๒๓๙-จ-๗๒๓๔
๑๘) นางสาวกนิษฐา เจริญเชื้อ	ทะเบียนเลขที่	ว-๒๓๙-จ-๗๒๓๖
๑๙) นายจิรกร สิมศิลา	ทะเบียนเลขที่	ว-๒๓๙-จ-๗๒๓๗
๒๐) นายชนาธิป สิงห์เกษมศักดิ์	ทะเบียนเลขที่	ว-๒๓๙-จ-๗๒๓๘
๒๑) นายวัชรกานต์ ประมาคยะเด	ทะเบียนเลขที่	ว-๒๓๙-จ-๗๒๔๐
๒๒) นายชอง เองขวัญกุล	ทะเบียนเลขที่	ว-๒๓๙-จ-๗๒๔๒
๒๓) นางสาวกฤษณา จันทุม	ทะเบียนเลขที่	ว-๒๓๙-จ-๗๒๔๒
๒๔) นางสาวพรนภา บุตรธรรม	ทะเบียนเลขที่	ว-๒๓๙-จ-๗๒๔๓
๒๕) นางสาวธาริณี อาจปลิว	ทะเบียนเลขที่	ว-๒๓๙-จ-๗๒๔๔
๒๖) นายธนโชติ ช่างล้อ	ทะเบียนเลขที่	ว-๒๓๙-จ-๗๒๔๖
๒๗) นางสาวพัชรา สมานันท์	ทะเบียนเลขที่	ว-๒๓๙-จ-๘๑๘๓
๒๘) นางสาวจุฑาทิพย์ แจ่มเรือน	ทะเบียนเลขที่	ว-๒๓๙-จ-๘๑๘๓
๒๙) นางสาวจณิสตา กุ้ยอ่อน	ทะเบียนเลขที่	ว-๒๓๙-จ-๘๑๘๗
๓๐) นางสาววรัญญา เขียนมัน	ทะเบียนเลขที่	ว-๒๓๙-จ-๘๑๘๘
๓๑) นางสาวจิรารัตน์ นุริตมบต์	ทะเบียนเลขที่	ว-๒๓๙-จ-๘๑๘๙

เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท ซีคอฟ จำกัด

เลขทะเบียน ว-๒๓๙

ที่ อก ๐๓๓๐(๑)/ ๕ ๑ ๘ ๐ ๔

ลงวันที่ ๒๑ ตุลาคม ๒๕๖๓

ขอขยายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๓๕๒ รายการ

น้ำเสีย จำนวน 46 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Aldrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
2	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[4] 2) Digestion, Inductively Coupled Plasma Method ^[4]
3	Barium	1) Digestion, Direct Nitrous Oxide-Acetylene Flame Method ^[4] 2) Digestion, Inductively Coupled Plasma Method ^[4]
4	α-BHC	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
5	β-BHC	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
6	γ-BHC	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
7	δ-BHC	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
8	Biochemical Oxygen Demand	1) 5-Day BOD Test, Azide Modification Method ^[4] 2) 5-Day BOD Test, Membrane Electrode Method ^[4]
9	Cadmium	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Method ^[4]

วิทย์

10 Chemical...

(นางริกาญจน์ ด้ตรสกลวิไล)

ผู้อำนวยการกลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
และทะเบียนห้องปฏิบัติการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
10	Chemical Oxygen Demand	1) Open Reflux, Titrimetric method ^[4] 2) Close Reflux, Colorimetric method ^[4] 3) Closed Reflux, Titrimetric Method ^[4]
11	Chlordane	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
12	Chromium	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Method ^[4]
13	Color	ADMI Weighted-Ordinate Spectrophotometric Method ^[4]
14	Copper	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Method ^[4]
15	Cyanide	Distillation, Colorimetric method ^[4]
16	2,4-D	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
17	4,4'-DDD	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
18	4,4'-DDE	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
19	4,4'-DDT	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
20	Dieldrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]

วิมล

21 Endosulfan I...

(นางริกาญจน์ ฉัตรสกุลวิไล)

ผู้อำนวยการกลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
และทะเบียนห้องปฏิบัติการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
21	Endosulfan I	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
22	Endosulfan II	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
23	Endosulfan Sulfate	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
24	Endrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
25	Endrin Aldehyde	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
26	Formaldehyde	Distillation, Colorimetric Method ^[3]
27	Free Chlorine	1) Iodometric Method ^[4] 2) DPD Colorimetric Method ^[4]
28	Heptachlor	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
29	Heptachlor epoxide	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
30	Hexavalent Chromium	1) Colorimetric Method ^[4] 2) Extraction, Air-Acetylene Flame Method ^[4]
31	Lead	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Method ^[4]

วิมล

(นางริกาญจน์ ฉัตรสกุลวิไล)

ผู้อำนวยการกลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ

32 Manganese...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
32	Manganese	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Method ^[4]
33	Mercury	Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[4]
34	Methoxychlor	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
35	Nickel	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Method ^[4]
36	Oil & Grease	1) Liquid-Liquid, Partition-Gravimetric Method ^[4] 2) Soxhlet Extraction Method ^[4]
37	pH	Electrometric Method ^[4]
38	Phenols	1) Distillation, Chloroform Extraction Method ^[4] 2) Distillation, Direct Photometric Method ^[4]
39	Selenium	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[4] 2) Digestion, Inductively Coupled Plasma Method ^[4]
40	Sulfide	1) Iodometric method ^[4] 2) Methylene blue method ^[4]
41	Temperature	Laboratory and Field Methods ^[4]
42	Total Dissolved Solids	Dried at 180 °C ^[4]
43	Total Kjeldahl Nitrogen	1) Macro Kjeldahl Method ^[4] 2) Semi-Micro Kjeldahl Method ^[4]
44	Total Suspended Solids	Dried at 103-105 °C ^[4]
45	Trivalent Chromium	1) Digestion, Direct Air-Acetylene Flame Method; Colorimetric Method; Calculation ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method; Colorimetric Method; Calculation ^[4] 3) Digestion, Inductively Coupled Plasma Method; Colorimetric Method; Calculation ^[4]
46	Zinc	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Method ^[4]

วิมล

(นางวิภาณูจน์ ฉัตรสกุลวิไล)

ผู้อำนวยการกลุ่มมาตรฐานวิธีการวิเคราะห์มลพิษ

น้ำใต้ดิน...

น้ำใต้ดิน จำนวน 123 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Acenaphthene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
2	Acetone	Purge and Trap Gas Chromatographic/ Mass Spectrometric Method ^[4]
3	Aldrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
4	Anthracene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
5	Antimony	Digestion, Inductively Coupled Plasma Spectrometric Method ^[4]
6	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[4] 2) Digestion, Inductively Coupled Plasma Method ^[4]
7	Atrazine	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
8	Barium	1) Digestion, Direct Nitrous Oxide-Acetylene Flame Method ^[4] 2) Digestion, Inductively Coupled Plasma Spectrometric Method ^[4]
9	Benz(a)anthracene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
10	Benzene	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]
11	Benzo(b)fluoranthene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
12	Benzo(k)fluoranthene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
13	Benzoic acid	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
14	Benzo(a)pyrene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
15	Benzo(g,h,i)perylene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]

วิมล

(นางวิภาณูจน์ ฉัตรสกุลวิไล)

ผู้อำนวยการกลุ่มมาตรฐานวิธีการวิเคราะห์มลพิษ

16 Beryllium...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
16	Beryllium	Digestion, Inductively Coupled Plasma Spectrometric Method ^[4]
17	Bis(2-chloroethyl)ether	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
18	Bis(2-ethylhexyl)phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
19	Bromodichloromethane	Purge and Trap Gas Chromatographic/Mass spectrometric Method ^[4]
20	Bromoform	Purge and Trap Gas Chromatographic/Mass spectrometric Method ^[4]
21	Butanol	Purge and Trap Gas Chromatographic/Mass spectrometric Method ^[4]
22	Butyl benzyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
23	Cadmium	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Spectrometric Method ^[4]
24	Carbazole	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
25	Carbon disulfide	Purge and Trap Gas Chromatographic/Mass spectrometric Method ^[4]
26	Carbon tetrachloride	Purge and Trap Gas Chromatographic/Mass spectrometric Method ^[4]
27	Chlordane	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
28	p-Chloroaniline	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
29	Chlorobenzene	Purge and Trap Gas Chromatographic/Mass spectrometric Method ^[4]
30	Chlorodibromomethane	Purge and Trap Gas Chromatographic/Mass spectrometric Method ^[4]
31	Chloroform	Purge and Trap Gas Chromatographic/Mass spectrometric Method ^[4]

(นางริภาณูจน์ ฉัตรสกุลวิไล)

ผู้อำนวยการศูนย์มาตรฐานวิธีการวิเคราะห์มลพิษ

กรมควบคุมมลพิษ

32 2-Chlorophenol...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
32	2-Chlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
33	Chromium	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Spectrometric Method ^[4]
34	Chromium (III)	1) Digestion, Direct Air-Acetylene Flame Method; Colorimetric Method; Calculation ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method; Colorimetric Method; Calculation ^[4] 3) Digestion, Inductively Coupled Plasma Spectrometric Method; Colorimetric Method; Calculation ^[4]
35	Chromium (VI)	1) Colorimetric Method ^[4] 2) Extraction, Air-Acetylene Flame Method ^[4]
36	Chrysene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
37	Cyanide	1) Distillation, Titrimetric Method ^[4] 2) Distillation, Colorimetric Method ^[4]
38	2,4-D	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
39	DDD	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
40	DDE	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
41	DDT	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]

(นางริภาณูจน์ ฉัตรสกุลวิไล)

ผู้อำนวยการศูนย์มาตรฐานวิธีการวิเคราะห์มลพิษ
และทะเบียนห้องปฏิบัติการ

42 Dibenz(a,h)...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
42	Dibenz(a,h)anthracene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
43	Di-n-butyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
44	1,2-Dichlorobenzene	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ⁽⁴⁾
45	1,3-Dichlorobenzene	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ⁽⁴⁾
46	1,4-Dichlorobenzene	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ⁽⁴⁾
47	3,3'-Dichlorobenzidine	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
48	1,1-Dichloroethane	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ⁽⁴⁾
49	1,2-Dichloroethane	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ⁽⁴⁾
50	1,1-Dichloroethylene	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ⁽⁴⁾
51	cis-1,2-Dichloroethylene	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ⁽⁴⁾
52	trans-1,2-Dichloroethylene	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ⁽⁴⁾
53	2,4-Dichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
54	1,2-Dichloropropane	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ⁽⁴⁾
55	1,3-Dichloropropane	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ⁽⁴⁾
56	1,3-Dichloropropene	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ⁽⁴⁾
57	Dieldrin	1) Liquid-Liquid Extraction, Gas Chromatographic/ Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
58	Diethyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾

อีก

59 2,4-Dimethylphenol...

(นางริกาญจน์ นิตกรกุลวิไล)

ผู้อำนวยการกลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
และทะเบียนห้องปฏิบัติการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
59	2,4-Dimethylphenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
60	2,4-Dinitrophenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
61	2,4-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
62	2,6-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
63	Di-n-Octyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
64	Endosulfan	1) Liquid-Liquid Extraction, Gas Chromatographic/ Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
65	Endrin	1) Liquid-Liquid Extraction, Gas Chromatographic/ Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
66	Ethylbenzene	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ⁽⁴⁾
67	Fluoranthene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
68	Fluorene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
69	Heptachlor	1) Liquid-Liquid Extraction, Gas Chromatographic/ Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
70	Heptachlor epoxide	1) Liquid-Liquid Extraction, Gas Chromatographic/ Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
71	Hexachlorobenzene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
72	Hexachloro-1,3-butadiene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾

อีก

73 n-Hexane...

(นางริกาญจน์ นิตกรกุลวิไล)

ผู้อำนวยการกลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
และทะเบียนห้องปฏิบัติการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
73	n-Hexane	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]
74	α -HCH	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
75	β -HCH	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
76	γ -HCH	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
77	Hexachlorocyclopentadiene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
78	Hexachloroethane	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
79	Indeno(1,2,3-cd)pyrene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
80	Isophorone	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
81	Lead	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Spectrometric Method ^[4]
82	Manganese	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Spectrometric Method ^[4]
83	Mercury	Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[4]
84	Methanol	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]

วิมล

85 Methoxychlor...

(นางริกาญจน์ ฉัตรสกุลวิไล)

ผู้อำนวยการกลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
และประเมินสิ่งแวดล้อม

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
85	Methoxychlor	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
86	Methyl bromide	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]
87	Methylene chloride	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]
88	2-Methylphenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
89	2-Methylnaphthalene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
90	Methyl tert-butyl ether	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]
91	Naphthalene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
92	Nickel	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Spectrometric Method ^[4]
93	Nitrobenzene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
94	N-Nitrosodiphenylamine	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
95	Polychlorinated Biphenyls - PCB-1016 - PCB-1221 - PCB-1232 - PCB-1242 - PCB-1248 - PCB-1254 - PCB-1260	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
96	Pentachlorophenol	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]

วิมล

(นางริกาญจน์ ฉัตรสกุลวิไล)

ผู้อำนวยการกลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
และประเมินสิ่งแวดล้อม

97 pH...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
97	pH	Electrometric method ^[4]
98	Phenanthrene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
99	Phenol	1) Distillation, Chloroform Extraction Method ^[4] 2) Distillation, Direct Photometric Method ^[4] 3) Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
100	Pyrene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
101	Selenium	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[4] 2) Digestion, Inductively Coupled Plasma Method ^[4]
102	Silver	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Inductively Coupled Plasma Method ^[4]
103	Styrene	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]
104	1,1,2,2-Tetrachloroethane	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]
105	Tetrachloroethylene	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]
106	Toluene	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]
107	TPH (C ₅ -C ₆)	Purge and Trap, Gas Chromatographic/ Mass spectrometric Method ^[7,9]
108	TPH (C ₈ -C ₁₆)	1) Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[6,8] 2) Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass spectrometric Method ^[6,9]
109	TPH (C ₁₆ -C ₃₅)	1) Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[6,8] 2) Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass spectrometric Method ^[6,9]
110	1,2,4-Trichlorobenzene	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]
111	1,1,1-Trichloroethane	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]

วิมล

112 1,1,2-Trichloroethane...

(นางริกาญจน์ ด้ตรสกุสวโร)

ผู้อำนวยการกลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ

กรมควบคุมมลพิษ กรุงเทพมหานคร

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
112	1,1,2-Trichloroethane	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]
113	Trichloroethylene	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]
114	2,4,5-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
115	2,4,6-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
116	1,3,5-Trimethylbenzene	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]
117	Vanadium	Digestion, Inductively Coupled Plasma Spectrometric Method ^[4]
118	Vinyl chloride	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]
119	m-Xylene	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]
120	o-Xylene	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]
121	p-Xylene	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]
122	Xylene (Total)	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]
123	Zinc	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Spectrometric Method ^[4]

อากาศเสีย (ปล่องระบาย) จำนวน 27 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Antimony	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]

วิมล

2 Arsenic...

(นางริกาญจน์ ด้ตรสกุสวโร)

ผู้อำนวยการกลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
และทะเบียนข้อมูลปฏิบัติการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
2	Arsenic	1) Isokinetic Sampling, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
3	Beryllium	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
4	Cadmium	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
5	Carbon monoxide	Instrumental Analyzer Method ^[5]
6	Chlorine	1) Absorption Sampling, Ion Chromatographic Method ^[5] 2) Isokinetic Sampling, Ion Chromatographic Method ^[5]
7	Chromium	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
8	Cobalt	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
9	Copper	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
10	Cresol	Adsorption Sampling, Gas Chromatographic Method ^[5]
11	Dioxin/Furans	Isokinetic Sampling, Analysis by ISO/IEC 17025 Accredited Laboratory or Analysis by Department of Industrial Works Registered Laboratory (Dioxins/Furans Analysis Approved) ^[5]
12	Hydrogen chloride	1) Absorption Sampling, Ion Chromatographic Method ^[5] 2) Isokinetic Sampling, Ion Chromatographic Method ^[5]
13	Hydrogen Fluoride	1) Absorption Sampling, Ion Chromatographic Method ^[5] 2) Isokinetic Sampling, Ion Chromatographic Method ^[5]



(นางริกาญจน์ จิตรสกุลไชย)

ผู้อำนวยการศูนย์มาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
และทะเบียนห้องปฏิบัติการ

14 Hydrogen Sulfide...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
14	Hydrogen Sulfide	Absorption Sampling, Iodometric Method ^[5]
15	Lead	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
16	Manganese	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
17	Mercury	Isokinetic Sampling, Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[5]
18	Nickel	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
19	Opacity	Ringelmann's Method ^[2]
20	Oxide of Nitrogen	1) Absorption Sampling, Ion Chromatographic Method ^[5] 2) Absorption Sampling, Phenoldisulfonic acid Method ^[5] 3) Instrumental Analyzer Method ^[5]
21	Selenium	1) Isokinetic Sampling, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
22	Sulfur dioxide	1) Isokinetic Sampling, Barium-Thorin Titrimetric Method ^[5] 2) Instrumental Analyzer Method ^[5]
23	Sulfuric acid	Isokinetic Sampling, Barium-Thorin Titrimetric Method ^[5]
24	Tin	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
25	Total Suspended Particulate	Isokinetic Sampling, Gravimetric Method ^[5]



(นางริกาญจน์ จิตรสกุลไชย)

ผู้อำนวยการศูนย์มาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
และทะเบียนห้องปฏิบัติการ

26 Vanadium...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
26	Vanadium	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
27	Xylene	1) Adsorption Sampling, Gas Chromatographic Method ^[6] 2) Adsorption Sampling, Gas Chromatographic/Mass Spectrometric Method ^[5]

สิ่งปฏิกูลหรือวัสดุที่ไม่ใช้แล้ว จำนวน 34 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Aldrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,9,22] 2) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 3) Soxhlet Extraction, Gas Chromatographic Method ^[10,22] 4) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26]
2	Antimony	1) Waste Extraction, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 3) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[7,16] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]
3	Arsenic	1) Waste Extraction, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 3) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[7,16] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]
4	Barium	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^[1,6,15] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14]

วิมล

(นางริกาญจน์ ฉัตรสกุลวิไล)

ผู้อำนวยการกลุ่มมาตรฐานวิธีการตรวจวัดมลพิษ
และประเมินผลกระทบต่อสุขภาพ

3) Digestion...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
5	Beryllium	3) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 4) Digestion, Inductively Coupled Plasma Method ^[7,14] 1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
6	Cadmium	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^[1,6,15] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 3) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]
7	Chlordane	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,9,22] 2) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 3) Soxhlet Extraction, Gas Chromatographic Method ^[10,22] 4) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26]
8	Chromium	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^[1,6,15] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 3) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]
9	Chromium (III)	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method; Waste Extraction, Colorimetric Method; Calculation Method ^[1,6,15,17] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method; Waste Extraction, Colorimetric Method; Calculation Method ^[1,6,16,17]

วิมล

(นางริกาญจน์ ฉัตรสกุลวิไล)

ผู้อำนวยการกลุ่มมาตรฐานวิธีการตรวจวัดมลพิษ
และประเมินผลกระทบต่อสุขภาพ

3) Digestion...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
10	Chromium (VI)	3) Digestion, Flame Atomic Absorption Spectrometric Method; Alkaline Digestion, Colorimetric Method; Calculation Method ^[7,8,15,17] 4) Digestion, Inductively Coupled Plasma Method; Alkaline Digestion, Colorimetric Method; Calculation Method ^[7,8,14,17]
11	Cobalt	1) Waste Extraction, Colorimetric Method ^[1,17] 2) Alkaline Digestion, Colorimetric Method ^[8,17]
12	Copper	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
13	2,4-D	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^[1,6,15] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 3) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]
14	DDD	1) Waste Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,24] 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[24]
15	DDE	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,9,22] 2) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 3) Soxhlet Extraction, Gas Chromatographic Method ^[10,22] 4) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26]



(นางวิภาณูจน์ ชิตรสกุลไธ)

ผู้อำนวยการศูนย์มาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
และทะเบียนห้องปฏิบัติการ

3) Soxhlet...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
16	DDT	3) Soxhlet Extraction, Gas Chromatographic Method ^[10,22] 4) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26] 1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,9,22] 2) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26]
17	Dieldrin	3) Soxhlet Extraction, Gas Chromatographic Method ^[10,22] 4) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26] 1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,9,22] 2) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26]
18	Endrin	3) Soxhlet Extraction, Gas Chromatographic Method ^[10,22] 4) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26] 1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,9,22] 2) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26]
19	Heptachlor	3) Soxhlet Extraction, Gas Chromatographic Method ^[10,22] 4) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26] 1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,9,22] 2) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 3) Soxhlet Extraction, Gas Chromatographic Method ^[10,22]



(นางวิภาณูจน์ ชิตรสกุลไธ)

ผู้อำนวยการศูนย์มาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
และทะเบียนห้องปฏิบัติการ

4) Soxhlet...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
20	Lead	4) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[9,26] 1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^[1,6,15] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 3) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]
21	Lindane	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,9,22] 2) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 3) Soxhlet Extraction, Gas Chromatographic Method ^[10,22] 4) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26]
22	Mercury	1) Waste Extraction, Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[1,18] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 3) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[19] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]
23	Methoxychlor	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,9,22] 2) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 3) Soxhlet Extraction, Gas Chromatographic Method ^[10,22] 4) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26]
24	Molybdenum	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]

วิมล

25 Nickel...

(นางริกาญจน์ ฉัตรสกุลวิไล)

ผู้อำนวยการกลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
และทะเบียนห้องปฏิบัติการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
25	Nickel	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^[1,6,15] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 3) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]
26	Polychlorinated Biphenyls - Aroclor 1016 - Aroclor 1221 - Aroclor 1232 - Aroclor 1242 - Aroclor 1248 - Aroclor 1254 - Aroclor 1260	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,9,23] 2) Soxhlet Extraction, Gas Chromatographic Method ^[10,23]
27	Pentachlorophenol	1) Waste Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[1,24] 2) Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[24]
28	pH	Electrometric Method ^[30,31]
29	Selenium	1) Waste Extraction, Digestion, Hydride Generation/ Atomic Absorption Spectrometric Method ^[1,6,20] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 3) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[7,20] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]
30	Silver	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
31	Thallium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
32	Trichloroethylene	1) Waste Extraction, Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[1,12,25] 2) Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[12,25]

วิมล

33 Vanadium...

(นางริกาญจน์ ฉัตรสกุลวิไล)

ผู้อำนวยการกลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
และทะเบียนห้องปฏิบัติการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
33	Vanadium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14]
34	Zinc	2) Digestion, Inductively Coupled Plasma Method ^[7,14] 1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^[1,6,15] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 3) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]

ดิน จำนวน 122 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Acenaphthene	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26]
2	Acetone	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
3	Aldrin	1) Ultrasonic Extraction, Gas Chromatographic Method ^[11,22] 2) Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
4	Anthracene	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26]
5	Antimony	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[7,15] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
6	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[7,15] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
7	Atrazine	Ultrasonic Extraction, Gas Chromatographic Method ^[11,22]
8	Barium	1) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]

9 Benz(a)anthracene...

(นางรวิภาญจน์ ด้ตรสกุลวิไล)

ผู้อำนวยการศูนย์มาตรฐานวิธีการวิเคราะห์ภาคตอนบนมลพิษ
และทะเบียนข้อมูลปฏิบัติการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
27	Chlordane	1) Ultrasonic Extraction, Gas Chromatographic Method ^[11,22] 2) Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
28	p-Chloroaniline	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26]
29	Chlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
30	Chlorodibromomethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
31	Chloroform	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
32	2-Chlorophenol	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26]
33	Chromium	1) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
34	Chromium (III)	1) Digestion, Flame Atomic Absorption Spectrometric Method; Colorimetric Method; Calculation Method ^[7,8,15,17] 2) Digestion, Inductively Coupled Plasma Method; Colorimetric Method; Calculation Method ^[7,8,14,17]
35	Chromium (VI)	Alkaline Digestion, Colorimetric Method ^[8,17]
36	Chrysene	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26]
37	Cyanide	1) Extraction, Distillation, Titrimetric Method ^[27,28,29] 2) Extraction, Distillation, Colorimetric Method ^[27,28,29]
38	2,4-D	Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[24]
39	DDD	1) Ultrasonic Extraction, Gas Chromatographic Method ^[11,22] 2) Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
40	DDE	1) Ultrasonic Extraction, Gas Chromatographic Method ^[11,22] 2) Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]

(นางรวิภาญจน์ ด้ตรสกุลวิไล)

ผู้อำนวยการศูนย์มาตรฐานวิธีการวิเคราะห์ภาคตอนบนมลพิษ
และทะเบียนข้อมูลปฏิบัติการ

41 DDT...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
41	DDT	1) Ultrasonic Extraction, Gas Chromatographic Method ^[11,22] 2) Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
42	Dibenz(a,h)anthracene	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26]
43	Di-n-butyl phthalate	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26]
44	1,2-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
45	1,3-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
46	1,4-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
47	3,3'-Dichlorobenzidine	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26]
48	1,1-Dichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
49	1,2-Dichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
50	1,1-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
51	cis-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
52	trans-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
53	2,4-Dichlorophenol	Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
54	1,2-Dichloropropane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
55	1,3-Dichloropropane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
56	1,3-Dichloropropene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]



57 Dieldrin...

(นางริกาญจน์ ฉัตรสกุลไชย)

ผู้อำนวยการศูนย์มาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
และทะเบียนห้องปฏิบัติการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
57	Dieldrin	1) Ultrasonic Extraction, Gas Chromatographic Method ^[11,22] 2) Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
58	Diethyl phthalate	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26]
59	2,4-Dimethylphenol	Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
60	2,4-Dinitrophenol	Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
61	2,4-Dinitrotoluene	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26]
62	2,6-Dinitrotoluene	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26]
63	Di-n-Octyl phthalate	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26]
64	Endosulfan	1) Ultrasonic Extraction, Gas Chromatographic Method ^[11,22] 2) Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
65	Endrin	1) Ultrasonic Extraction, Gas Chromatographic Method ^[11,22] 2) Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
66	Ethylbenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
67	Fluoranthene	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26]
68	Fluorene	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26]
69	Heptachlor	1) Ultrasonic Extraction, Gas Chromatographic Method ^[11,22] 2) Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]



70 Heptachlor epoxide...

(นางริกาญจน์ ฉัตรสกุลไชย)

ผู้อำนวยการศูนย์มาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
และทะเบียนห้องปฏิบัติการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
70	Heptachlor epoxide	1) Ultrasonic Extraction, Gas Chromatographic Method ^[11,22] 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26]
71	Hexachlorobenzene	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26]
72	Hexachloro-1,3-butadiene	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26]
73	n-Hexane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,25]
74	α -HCH	1) Ultrasonic Extraction, Gas Chromatographic Method ^[11,22] 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26]
75	β -HCH	1) Ultrasonic Extraction, Gas Chromatographic Method ^[11,22] 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26]
76	γ -HCH	1) Ultrasonic Extraction, Gas Chromatographic Method ^[11,22] 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26]
77	Hexachlorocyclopentadiene	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26]
78	Hexachloroethane	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26]
79	Indeno(1,2,3-cd)pyrene	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26]
80	Isophorone	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26]
81	Lead	1) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
82	Manganese	1) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]



(นางริกาญจน์ ฉัตรสกุลไชย)

ผู้อำนวยการศูนย์มาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
และทะเบียนห้องปฏิบัติการ

83 Mercury...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
83	Mercury	1) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[19] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
84	Methanol	Ultrasonic Extraction, Direct Aqueous Injection, Gas Chromatographic Method ^[11,21]
85	Methoxychlor	1) Ultrasonic Extraction, Gas Chromatographic Method ^[11,22] 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26]
86	Methyl bromide	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,25]
87	Methylene chloride	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,25]
88	2-Methylphenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26]
89	2-Methylnaphthalene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26]
90	Methyl tert-butyl ether	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,25]
91	Naphthalene	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26]
92	Nickel	1) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
93	Nitrobenzene	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26]
94	N-Nitrosodiphenylamine	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26]
95	Polychlorinated Biphenyls - Aroclor 1016 - Aroclor 1221 - Aroclor 1232 - Aroclor 1242 - Aroclor 1248 - Aroclor 1254 - Aroclor 1260	Soxhlet Extraction, Gas Chromatographic Method ^[10,23]



(นางริกาญจน์ ฉัตรสกุลไชย)

ผู้อำนวยการศูนย์มาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
และทะเบียนห้องปฏิบัติการ

96 Pentachlorophenol...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
96	Pentachlorophenol	Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[24]
97	Phenanthrene	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26]
98	Phenol	Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
99	Pyrene	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26]
100	Selenium	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[7,20] 2) Digestion, Inductively Coupled Plasma Method ^[7,16]
101	Silver	1) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
102	Styrene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
103	1,1,2,2-Tetrachloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
104	Tetrachloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
105	Toluene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
106	TPH (C ₅ -C ₈)	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
107	TPH (C ₈ -C ₁₆)	1) Soxhlet Extraction, Gas Chromatographic Method ^[10,21] 2) Soxhlet Extraction, Gas Chromatographic/ Mass spectrometric Method ^[10,21]
108	TPH (C ₁₆ -C ₃₅)	1) Soxhlet Extraction, Gas Chromatographic Method ^[10,21] 2) Soxhlet Extraction, Gas Chromatographic/ Mass spectrometric Method ^[10,23]
109	1,2,4-Trichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
110	1,1,1-Trichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]

111 1,1,2-Trichloroethane...

(นางริกาญจน์ อัครสกุลวิไล)

ผู้อำนวยการศูนย์มาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
และทะเบียนต้องปฏิบัติการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
111	1,1,2-Trichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
112	Trichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
113	2,4,5-Trichlorophenol	Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
114	2,4,6-Trichlorophenol	Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
115	1,3,5-Trimethylbenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
116	Vanadium	Digestion, Inductively Coupled Plasma Method ^[7,14]
117	Vinyl chloride	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
118	m-Xylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
119	o-Xylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
120	p-Xylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
121	Xylene (Total)	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
122	Zinc	1) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]

เอกสารอ้างอิง

- กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2548. เรื่อง การกำจัดสิ่งปฏิกูลหรือ
วัสดุที่ไม่ใช้แล้ว. ราชกิจจานุเบกษา. 25 มกราคม 2549. เล่มที่ 123 ตอนพิเศษ 11ง.
- กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2549. เรื่อง กำหนดค่าปริมาณเขม่า
ควันที่เจือปนในอากาศที่ระบายออกจากปล่องของหม้อน้ำโรงสีข้าวที่ใช้กลบเป็นเชื้อเพลิง.
ราชกิจจานุเบกษา. 4 ธันวาคม 2549. เล่มที่ 123 ตอนพิเศษ 125ง.
- สมาคมวิศวกรรมสิ่งแวดล้อมแห่งประเทศไทย. คู่มือวิเคราะห์น้ำเสีย. พิมพ์ครั้งที่ 4. กรุงเทพฯ:
เรือนแก้วการพิมพ์, 2547.
- APHA, AWWA, WEF. Standard Methods for the Examination of Water and
Wastewater. 23rd ed. Washington, DC: APHA, 2017.
- United States Environmental Protection Agency. Standards of Performance for
New Stationary Sources. 40 CFR 60. Appendix A, 2019.

6. United States...

(นางริกาญจน์ อัครสกุลวิไล)

ผู้อำนวยการศูนย์มาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
และทะเบียนต้องปฏิบัติการ

6. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. SW-846, 1997.
7. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Acid Digestion of Sediments, Sludges, and Soils. SW-846 Method 3050B, 1996.
8. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Alkaline Digestion for Hexavalent Chromium. SW-846 Method 3060A, 1996.
9. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Separatory Funnel Liquid-Liquid Extraction. SW-846 Method 3510C, 1996.
10. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Soxhlet Extraction. SW-846 Method 3540C, 1996.
11. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Ultrasonic Extraction. SW-846 Method 3550C, 2007.
12. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Purge-and-Trap for Aqueous Samples. SW-846 Method 5030C, 2003.
13. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Closed-System Purge-and-Trap And Extraction For Volatile Organics in Soil and Waste Samples. SW-846 Method 5035A, 2002.
14. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Inductively Coupled Plasma-optical Emission Spectrometry. SW-846 Method 6010D, 2018
15. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Flame Atomic Absorption Spectrophotometry. SW-846 Method 7000B, 2007.
16. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Antimony and Arsenic (Atomic Absorption, Borohydride Reduction). SW-846 Method 7062, 1992.
17. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Chromium, Hexavalent (Colorimetric), SW-846 Method 7196A, 1992.
18. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Mercury in Liquid Waste (Manual Cold-Vapor Technique, SW-846 Method 7470A, 1994.
19. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Mercury in Solid or Semisolid Waste (Manual Cold-Vapor Technique, SW-846 Method 7471B, 2007.

(นางริการัญณ์ ฉัตรสกุลวิไล)

ผู้อำนวยการกลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
กรมควบคุมมลพิษ

20. United States...

20. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Selenium (Atomic Absorption, Borohydride Reduction), SW-846 Method 7742, 1994.
21. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Nonhalogenated Organics Using GC/FID. SW-846 Method 8015D, 2003.
22. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Organochlorine Pesticide by Gas Chromatography. SW-846 Method 8081B, 2007.
23. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Polychlorinated Biphenyls (PCBs) By Gas Chromatography. SW-846 Method 8082A, 2007.
24. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Chlorinated Herbicides By GC Using Methylation or Pentafluorobenzoylation Derivatization. SW-846 Method 8151A, 1996.
25. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Volatile Organic Compounds by Gas Chromatography/ Mass Spectrometry (GC/MS). SW-846 Method 8260D, 2018.
26. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. SemiVolatile Organic Compounds by Gas Chromatography/Mass Spectrometry. SW-846 Method 8270E, 2018.
27. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Total and Amenable Cyanide: Distillation. SW-846 Method 9010C, 2004.
28. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Cyanide Extraction Procedure for Solids and Oils. SW-846 Method 9013A, 2014.
29. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Cyanide in Waters and Extracts Using Titrimetric and Manual Spectrophotometric. SW-846 Method 9014, 2014.
30. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. pH Electrometric Measurement. SW-846 Method 9040C, 2004.
31. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Solid and Waste pH. SW-846 Method 9045D, 2004.

วิมล

(นางริการัญณ์ ฉัตรสกุลวิไล)

ผู้อำนวยการกลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษ
และทะเบียนห้องปฏิบัติการ

ภาคผนวก ซ

ใบรับรองความสามารถห้องปฏิบัติการและขอบข่าย
การรับรองห้องปฏิบัติการทดสอบ ตาม ISO/IEC 17025
จากสำนักงานมาตรฐานอุตสาหกรรม (สมอ.)



ใบรับรองเลขที่ 20T173/1151

ใบรับรองห้องปฏิบัติการ

อาศัยอำนาจตามความในพระราชบัญญัติการมาตรฐานแห่งชาติ พ.ศ. ๒๕๕๑

เลขานุการสำนักงานมาตรฐานผลิตภัณฑ์อุตสาหกรรม

ออกใบรับรองฉบับนี้ให้

บริษัท ซีคอฟ จำกัด

มีห้องปฏิบัติการตั้งอยู่เลขที่

239 ถนนริมคลองประปา แขวงบางซื่อ เขตบางซื่อ กรุงเทพมหานคร

ได้รับการรับรองความสามารถห้องปฏิบัติการทดสอบ

ตามมาตรฐานเลขที่ มอก. 17025-2561 (ISO/IEC 17025 : 2017)

ข้อกำหนดทั่วไปว่าด้วยความสามารถห้องปฏิบัติการทดสอบและสอบเทียบ

หมายเลขการรับรองที่ ทดสอบ ๐๓๙๔

โดยมีสาขาการรับรองตามรายละเอียดแนบท้ายใบรับรอง

ตั้งแต่วันที่ ๙ กันยายน พ.ศ. ๒๕๖๓

ถึง วันที่ ๘ กันยายน พ.ศ. ๒๕๖๖

ออกให้ ณ วันที่ ๒๓ กันยายน ๒๕๖๓

(นายระกิด์ รินกิจธนวิทย์)

รองเลขาธิการ ปฏิบัติราชการแทน

เลขาธิการสำนักงานมาตรฐานผลิตภัณฑ์อุตสาหกรรม



กระทรวงอุตสาหกรรม สำนักงานมาตรฐานผลิตภัณฑ์อุตสาหกรรม

รายละเอียดแนบท้ายใบรับรองห้องปฏิบัติการทดสอบ

ใบรับรองเลขที่ 20T173/1151

ชื่อห้องปฏิบัติการ

ที่อยู่

หมายเลขการรับรองที่

สถานที่ห้องปฏิบัติการ

ห้องปฏิบัติการทดสอบ บริษัท ซีคอฟ จำกัด

239 ถนนริมคลองประปา แขวงบางซื่อ เขตบางซื่อ กรุงเทพมหานคร

ทดสอบ 0394

☒ ภายใน ☐ นอกสถานที่ ☐ชั่วคราว ☐เคลื่อนที่

สาขาการทดสอบ	รายการทดสอบ	วิธีทดสอบ
สาขาสิ่งแวดล้อม		
1. น้ำและน้ำเสีย (water and wastewater)	- Arsenic 0.000 5 mg/l to 0.090 0 mg/l	- Standard Methods for the Examination of Water and Wastewater, APHA, AWWA, WEF, 23 rd edition, 2017, Part 3030 F and Part 3114 C
	- Arsenic 0.05 mg/l to 4.50 mg/l	- Standard Methods for the Examination of Water and Wastewater, APHA, AWWA, WEF, 23 rd edition, 2017, Part 3030 E and Part 3120 B
	- Barium 0.02 mg/l to 4.50 mg/l	
	- Cadmium 0.01 mg/l to 4.50 mg/l	
	- Chromium 0.01 mg/l to 4.50 mg/l	
	- Copper 0.02 mg/l to 4.50 mg/l	
	- Iron 0.05 mg/l to 9.00 mg/l	
	- Lead 0.03 mg/l to 4.50 mg/l	
	- Manganese 0.01 mg/l to 9.00 mg/l	
	- Nickel 0.01 mg/l to 4.50 mg/l	
	- Zinc 0.02 mg/l to 9.00 mg/l	

ฉบับที่ 1 ตั้งแต่วันที่ 9 กันยายน 2563

หน้า 1/5

กระทรวงอุตสาหกรรม สำนักงานมาตรฐานผลิตภัณฑ์อุตสาหกรรม

รายละเอียดแนบท้ายใบรับรองห้องปฏิบัติการทดสอบ

ใบรับรองเลขที่ 20T173/1151

หมายเลขการรับรองที่ ทดสอบ 0394

สถานภาพห้องปฏิบัติการ ☒ ถาวร ☐ นอกสถานที่ ☐ชั่วคราว ☐เคลื่อนที่

สาขาการทดสอบ	รายการทดสอบ	วิธีทดสอบ
<p>สาขาสังแวดล้อม</p> <p>1. น้ำและน้ำเสีย (ต่อ) (water and wastewater) (cont.)</p> <p>2. คุณภาพอากาศ (air quality)</p> <p>2.1 บริเวณทำงาน (workplace)</p>	<p>- COD 100 mg/l to 4 000 mg/l</p> <p>- Total dust 0.10 mg/filter to 2.00 mg/filter</p> <p>- Respirable dust 0.10 mg/filter to 2.00 mg/filter</p> <p>- Benzene 1.10 µg/tube to 420 µg/tube</p> <p>- Toluene 1.10 µg/tube to 420 µg/tube</p> <p>- Total xylenes 2.20 µg/tube to 840 µg/tube</p> <p>• m,p-xylene 1.10 µg/tube to 420 µg/tube</p> <p>• o-xylene 1.10 µg/tube to 420 µg/tube</p>	<p>- Standard Methods for the Examination of Water and Wastewater, APHA, AWWA, WEF, 23rd edition, 2017, Part 5220 D</p> <p>- NIOSH Manual of Analytical Methods (NMAM), method 0500, 4th edition, 15th August 1994 (Exclude Sampling)</p> <p>- NIOSH Manual of Analytical Method(NMAM), method 0600, 4th edition, 15th January 1998 (Exclude Sampling)</p> <p>- NIOSH Manual of Analytical Methods (NMAM) , method 1501, 4th edition, 15th March 2003 (Exclude Sampling)</p>

รายละเอียดแนบท้ายใบรับรองห้องปฏิบัติการทดสอบ

ใบรับรองเลขที่ 20T173/1151

หมายเลขการรับรองที่ ทดสอบ 0394

สถานภาพห้องปฏิบัติการ ☒ ถาวร ☐ นอกสถานที่ ☐ชั่วคราว ☐เคลื่อนที่

สาขาการทดสอบ	รายการทดสอบ	วิธีทดสอบ
<p>สาขาสังแวดล้อม</p> <p>2. คุณภาพอากาศ (ต่อ) (air quality) (cont.)</p> <p>2.2 อากาศในปล่องระบาย อากาศ (stack)</p> <p>2.3 บรรยากาศทั่วไป (ambient air)</p>	<p>- Sulfur dioxide 1.00 mg/l to 16 000 mg/l (solution)</p> <p>- Hydrogen fluoride 5 µg/sample to 400 µg/sample</p> <p>- Hydrogen chloride 5 µg/sample to 400 µg/sample</p> <p>- Volatile organic compounds (VOCs)</p> <p>• Chloroethene 0.05 µg/m³ to 51.00 µg/m³</p> <p>• 1,3 - butadiene 0.04 µg/m³ to 44.00 µg/m³</p> <p>• Bromomethane 0.08 µg/m³ to 77.00 µg/m³</p> <p>• Acrolein 0.05 µg/m³ to 45.00 µg/m³</p> <p>• Acrylonitrile 0.04 µg/m³ to 43.00 µg/m³</p> <p>• Dichloromethane 0.14 µg/m³ to 69.00 µg/m³</p> <p>• Carbon disulfide 0.06 µg/m³ to 62.00 µg/m³</p> <p>• Trichloromethane 0.20 µg/m³ to 97.00 µg/m³</p>	<p>- US.EPA , Code of Federal Regulations, 40 CFR 60 appendix A, Method 6, July 2019 (Exclude Sampling)</p> <p>- In-house method : WI-7.2-1-22 based on US.EPA, Code of Federal Regulations, 40 CFR 60 appendix A Method 26, 2019 (Exclude Sampling)</p> <p>- In-house method :WI-7.2-1-24 based on US.EPA , Compendium Method TO - 15, EPA / 625 / R-96 / 010b, January 1999 (Include sampling)</p>

รายละเอียดแนบท้ายใบรับรองห้องปฏิบัติการทดสอบ
ใบรับรองเลขที่ 20T173/1151

หมายเลขการรับรองที่ ทดสอบ 0394
สถานภาพห้องปฏิบัติการ ☒ ถาวร ☐ นอกสถานที่ ☐ชั่วคราว ☐เคลื่อนที่

สาขาการทดสอบ	รายการทดสอบ	วิธีทดสอบ
<p>สาขาสิ่งแวดล้อม</p> <p>2. คุณภาพอากาศ (ต่อ) (air quality) (cont.)</p> <p>2.3 บรรยากาศทั่วไป (ต่อ) (ambient air) (cont.)</p>	<p>- Volatile organic compounds (VOCs) (cont.)</p> <ul style="list-style-type: none"> 1,2 - dichloroethane 0.08 µg/m³ to 80.00 µg/m³ Benzene 0.06 µg/m³ to 63.00 µg/m³ Carbon tetrachloride 0.25 µg/m³ to 125 µg/m³ Trichloroethylene 0.21 µg/m³ to 107 µg/m³ 1,2 - dichloropropane 0.18 µg/m³ to 92.00 µg/m³ Tetrachloroethylene 0.27 µg/m³ to 135 µg/m³ 1,2 - dibromoethane 0.31 µg/m³ to 153 µg/m³ 1,1,2,2 - tetrachloroethane 0.69 µg/m³ to 137 µg/m³ 	<p>- In-house method :WI-7.2-1-24 US.EPA , Compendium Method TO - 15, EPA / 625 / R-96 / 010b, January 1999 (Include sampling)</p>

ฉบับที่ 1 ตั้งแต่ วันที่ 9 กันยายน 2563 หน้า 4/5
กระทรวงอุตสาหกรรม สำนักงานมาตรฐานผลิตภัณฑ์อุตสาหกรรม

รายละเอียดแนบท้ายใบรับรองห้องปฏิบัติการทดสอบ
ใบรับรองเลขที่ 20T173/1151

หมายเลขการรับรองที่ ทดสอบ 0394
สถานภาพห้องปฏิบัติการ ☒ ถาวร ☐ นอกสถานที่ ☐ชั่วคราว ☐เคลื่อนที่

สาขาการทดสอบ	รายการทดสอบ	วิธีทดสอบ
<p>สาขาสิ่งแวดล้อม</p> <p>2. คุณภาพอากาศ (ต่อ) (air quality) (cont.)</p> <p>2.3 บรรยากาศทั่วไป (ต่อ) (ambient air) (cont.)</p>	<p>- Volatile organic compounds (VOCs) (cont.)</p> <ul style="list-style-type: none"> Benzyl chloride 0.52 µg/m³ to 103 µg/m³ 1,4 - dichlorobenzene 0.24 µg/m³ to 120 µg/m³ 	<p>- In-house method :WI-7.2-1-24 US.EPA , Compendium Method TO - 15, EPA / 625 / R-96 / 010b, January 1999 (Include sampling)</p>

ออกให้ ณ วันที่ 9 กันยายน 2563

(นายวีระกิตต์ รินทกิจธนวิชัย)
รองเลขาธิการ ปฏิบัติราชการแทน
เลขาธิการสำนักงานมาตรฐานผลิตภัณฑ์อุตสาหกรรม

ฉบับที่ 1 ตั้งแต่ วันที่ 9 กันยายน 2563 หน้า 5/5
กระทรวงอุตสาหกรรม สำนักงานมาตรฐานผลิตภัณฑ์อุตสาหกรรม