

ภาคผนวก จ

เอกสารสอบเทียบเครื่องมือ

---

รายการใบรับรองสอบเทียบ/ทวนสอบ เครื่องมือหลักประจำห้องปฏิบัติการวิเคราะห์ สำหรับวิเคราะห์คุณภาพสิ่งแวดล้อม

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*	Remark
เครื่องมือหลักประจำห้องปฏิบัติการวิเคราะห์คุณภาพน้ำ									
1	pH Meter	ความเป็นกรดเป็นด่าง (pH)	Mettler-Toledo	Seven Easy S20 / 1231155210	National Food Institute, Ministry of Industry, Thailand	2301846-001-01	24 Feb 23	23 Feb 24	-
2	pH Meter		Mettler-Toledo	Seven Easy S20 / 1230525212	National Food Institute, Ministry of Industry, Thailand	2302181-001-01	24 Mar 23	22 Mar 24	-
3	Analytical Balance (Readability 0.01 mg)	ของแข็งแขวนลอยทั้งหมด (SS) ของแข็งละลายน้ำ (TDS)	Mettler-Toledo	XSR205DU / C210685394	Technology Promotion Association (Thailand-Japan)	23MM113	26 Apr 23	24 Apr 24	-
4	Hot Air Oven		Memmert	UF55 / B216.1666	Technology Promotion Association (Thailand-Japan)	22TM1490	19 Oct 22	18 Oct 23	-
5	Analytical Balance (Readability 0.1 mg)	น้ำมัน และไขมัน (Fat, Oil & Grease)	Mettler-Toledo	XSR204 / C117635043	National Food Institute, Ministry of Industry, Thailand	2302827-001-01	10 May 23	8 May 24	-
6	BOD Incubator	บีโอดี(BOD)	Arco	UC4-1320 / (UAE.WAO.015/2561)	Technology Promotion Association (Thailand-Japan)	23TM249	15 Feb 23	14 Feb 24	-
7	BOD Incubator		Arco	UR-1320 / (UAE.WAO.018/2551)	Technology Promotion Association (Thailand-Japan)	23TM375	12 Apr 23	10 Apr 24	-
8	Digester Unit	ไนโตรเจนในรูปทีเคเอ็น (TKN)	FOSS TECATOR	2520auto / 91794469	National Food Institute, Ministry of Industry, Thailand	2302413-001-01	30 Mar 23	28 Mar 24	-
9	Distillation Unit (Kjeldahl Method)		FOSS TECATOR	KT8100/ 91889052	FOSS South East Asia	6623	25 Jul 22	24 Jul 23	-
10	Incubator	โคลิฟอร์มแบคทีเรีย (Coliform Bacteria) ฟีคอลลีฟอร์มแบคทีเรีย	Binder	KB400 / 20200000015535	Technology Promotion Association (Thailand-Japan)	23TM726	27 Apr 23	25 Apr 24	-
11	Incubator		Memmert	IPP 260 / V616.0066	Technology Promotion Association (Thailand-Japan)	23TM728	27 Apr 23	25 Apr 24	-

รายการใบรับรองสอบเทียบ/ทวนสอบ เครื่องมือหลักประจำห้องปฏิบัติการวิเคราะห์ สำหรับวิเคราะห์คุณภาพสิ่งแวดล้อม

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*	Remark
เครื่องมือหลักประจำห้องปฏิบัติการวิเคราะห์คุณภาพน้ำ									
12	Water Bath		Memmert	WNE 14 / L416.0606	Technology Promotion Association (Thailand-Japan)	<a href="#">23TM193</a>	15 Feb 23	14 Feb 24	-
13	Water Bath		Memmert	WNE 14 / L416.0612	Technology Promotion Association (Thailand-Japan)	<a href="#">23TM194</a>	15 Feb 23	14 Feb 24	-
14	Auto Clave		ALP	CL-40L / 808763	Technology Promotion Association (Thailand-Japan)	<a href="#">23TM763</a>	27 Apr 23	25 Apr 24	-
15	Auto Clave		ALP	CL-40L / 810010	DKSH (Thailand) Ltd.	<a href="#">C11230106</a>	9 Jun 23	7 Jun 24	-
16	Refrigerator		Sanyo	SBC-337KD(GYN) / 71100607	Technology Promotion Association (Thailand-Japan)	<a href="#">23TM376</a>	11 Apr 23	9 Apr 24	-
17	Refrigerator		Sanden	YPM11P / YPM11P-150801805	Technology Promotion Association (Thailand-Japan)	<a href="#">22TM1067</a>	11 Jul 22	10 Jul 23	-
18	Analytical Balance		OHAUS	PX623 / C236754745	DKSH (Thailand) Ltd.	<a href="#">C01223732</a>	9 Dec 22	8 Dec 23	-

Due Date of Calibration\* : Based on the annual calibration plan. At least 1 time per year.

## Calibration Certificate

Certificate No.: 2301846-001-01  
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.  
Address: 3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchack, Prakanong, Bangkok 10260

Page 1 of 5

Equipment: pH Meter  
Manufacturer: Mettler Toledo  
Model: SevenEasy TM S20 pH  
Serial No.: 1231155210  
ID No.: UAE.WAT.010/2553  
Order No.: 2301846  
Operation No.: 2301846-001  
Date of Receipt: 17 February 2023  
Date of Calibration: 24 February 2023

Calibrated by Mr.Worapob Sooktong  
Scientist  
Approved by (Mr.Nattapong Niyomchan)  
Specialist, Division of Calibration Laboratory  
Responsible for the Technical Management Team  
Date of Issue: 24 February 2023

The uncertainties are for a confidence probability of approximately 95%.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65

## Calibration Report

Certificate No.: 2301846-001-01  
Equipment: pH Meter  
Resolution: 0.01 pH ; 1 mV  
Manufacturer: Mettler Toledo  
Model: SevenEasy TM S20 pH  
Serial No.: 1231155210  
Type: Bench top  
ID No.: UAE.WAT.010/2553

Date of Calibration: 24 February 2023 Page 2 of 5

Location: Chemical Calibration Laboratory, National Food Institute  
Environment Condition: Ambient Temperature: ( 25.1 ± 1.5 ) °C Relative Humidity: ( 50 ± 5 ) %  
Condition of Equipment: Good Condition  
Condition of this Results of Calibration

1. Calibration Method: In house method : W-CC-002 based on direct measurement by using standard voltage calibrator and certified reference material (CRM)

2. Reference Standards / Certified Reference Material

Instruments	Serial / ID No.	Manufacturer	Certificate No.	Due Date
2.1 DC Voltage Calibrator	2709007	Fluke	2281959	17 June 2023
2.2 Digital Thermometer	2709007	Fluke	CC 650677-01	30 October 2023
2.3 Thermo-Hygro Meter	NPI.BTH 007115	PONPE 450	QR22-0886	26 April 2023
Certified Reference Material				
	Lot No.	Manufacturer	Ref. No.	Expiry Date
2.4 pH buffer 4.006 (Primary pH buffer Solution)	832806	CPAchem	PH216.L5	8 August 2024
2.5 pH buffer 6.865 (Primary pH buffer Solution)	832807	CPAchem	PH217.L5	8 August 2024
2.6 pH buffer 10.01 (Primary pH buffer Solution)	832809	CPAchem	PH220.L5	8 August 2023
2.7 pH buffer 7.00 (Standard pH buffer Solution)	832610	CPAchem	PH107.L5	8 August 2023

3. This certification is traceable to The International System of Unit (SI Unit)

3.1 Instruments No 2.1	through	NSC-TIS-TIS 17025 Laboratory Accreditation of Calibration No.0008
3.2 Instruments No 2.2	through	NSC-TIS-TIS 17025 Laboratory Accreditation of Calibration No.0061
3.3 Instruments No 2.3	through	NSC-TIS-TIS 17025 Laboratory Accreditation of Calibration No.0292
3.4 Certified Reference Material No. 2.4 to 2.6	traceable to	Primary measurement method- formed call using calibrated thermometer, barometer, and nanovoltmeter. The Standard Solution preparation and certified by CPAchem Ltd is accredited to ISO 17034 and ISO/IEC 17025
3.5 Certified Reference Material No.2.7	traceable to	SI M Ref N 18-27 Lot N 04.05.2021; SI M Ref N 18-28 Lot N 28.05.2021; SI M Ref N 18-27 Lot N 04.05.2021; SI M Ref N 18-28 Lot N 28.05.2021, the Standard Solution preparation and certified by CPAchem Ltd is accredited to ISO 17034 and ISO/IEC 17025

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

F-CS-012 Revision: 01 Date: 20-04-65

## Calibration Report

Certificate No.: 2301846-001-01  
Equipment: pH Meter  
Resolution: 0.01 pH ; 1 mV  
Manufacturer: Mettler Toledo  
Model: SevenEasy TM S20 pH  
Serial No.: 1231155210  
Type: Bench top  
ID No.: UAE.WAT.010/2553

Date of Calibration: 24 February 2023 Page 3 of 5

Calibration Results:  
1. Calibration of pH Meter ( Manual Temperature Compensation at 25 °C )

Nominal pH	DC Voltage Standard (mV)	Average Indicator Reading		Uncertainty ( ± mV )	Coverage Factor ( k )
		mV	pH		
0	414.120	414	0.00	0.58	2.50
2	295.814	296	2.00	0.58	2.50
4	177.464	178	4.00	0.58	2.50
6	59.180	59	6.00	0.58	2.50
7	0.000	0	7.00	0.58	2.50
8	-59.158	-59	8.00	0.58	2.50
10	-177.460	-177	10.00	0.58	2.50
12	-295.811	-296	12.00	0.58	2.50
14	-414.117	-414	14.00	0.58	2.50

2. Calibration of pH Meter with Electrode ( Manual Temperature Compensation at 25 °C )

Equipment: pH Electrode  
Type: Combined Electrode  
Manufacturer: Mettler Toledo  
Model: InLab Solids  
Serial No.: 9015311  
ID No.: N/A  
Performance of Electrode system ( Three-Point Calibration at pH 4, pH 7 and pH 10 )

Certified Value @25 °C (pH)	Average Indicator Reading		Relative Slope (%)	Uncertainty ( ± pH )	Coverage Factor ( k )
	pH	mV			
4.008	4.01	186	-	0.0071	2.00
6.985	6.99	19	97.98	0.0075	2.00
10.008	10.01	-180	97.29	0.0095	2.00
6.985	6.99	15	-	0.0092	2.00

F-CS-012 Revision: 01 Date: 20-04-65

## Calibration Report

Certificate No.: 2301846-001-01  
Equipment: Digital Thermometer with RTD  
Resolution: 0.1 °C  
Model: SevenEasy TM S20 pH  
Serial No.: 1231155210  
ID No.: UAE.WAT.010/2553  
Manufacturer: Mettler Toledo

Date of Calibration: 24 February 2023 Page 4 of 5

Location: Chemical Calibration Laboratory, National Food Institute  
Environment Condition: Ambient Temperature 25 °C ± 1 °C  
Relative Humidity 48 % ± 3 %

Condition of this results of Calibration:

- Calibration Method : - In house method: W-TE-025 by comparison with standard thermometer.  
- The Calibration is determined by comparing with a known temperature from a standard resistance thermometer.  
- The temperature scale in use at this laboratory is the International Temperature scale of 1990 ( ITS-90 ).

2. Reference Standard Instrument :

Instrument	Model	Serial No.	Certificate No.	Due Date	Through
HANDHELD THERMOMETER	1523	2118154	PSL-T 0673/65	07-Jun-23	TISTR
Platinum Resistance Thermometer (PRT)	5627A	877332			

Support Equipment : - Low Temperature Bath (Micro Bath), Model: 7103, S/N: A39538,AN65 A85181.

- This certificate is traceable to International System of Units (SI Units).
- This certificate was certified only for the instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.
- Condition of Calibrated Item : Good
- Result of Calibration : ☒ Without adjustment ☐ After adjustment

F-CS-012 Revision: 01 Date: 20-04-65

## Calibration Report

**Certificate No.:** 2301846-001-01  
**Equipment:** Digital Thermometer with RTD  
Resolution: 0.1 °C Model: SevenEasy TM S20 pH  
Serial No.: 1231155210 ID No.: UAE.WAT.010/2553  
Manufacturer: Mettler Toledo

**Date of Calibration:** 24 February 2023 Page 5 of 5

**Calibration point:** 15.0, 25.0 and 35.0 °C  
**Calibration result:**  
- The probe was immersed in liquid bath or dry bath to a minimum depth of 120 mm.  
- Description of probe, model : S/N :  
Dimension of probe : Diameter 9 mm., Length 120 mm.,  
Sheath material : Stainless Steel

UUC* Reading (°C)	Standard Temperature (°C)	Correction Value (°C)	Uncertainty ± (°C)
15.1	15.015	- 0.1	0.11
25.0	25.014	0.0	0.11
35.1	35.016	- 0.1	0.11

**Note**  
- UUC\* : Unit Under Calibration

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k= 2, providing a level of confidence of approximately 95 %.

\*\*\*\*\* End \*\*\*\*\*

F-CS-012 Revision: 01 Date: 20-04-65

เอกสารไม่ควบคุม



## Calibration Certificate

**Certificate No.:** 2302181-001-01  
**Client name:** UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.  
**Address:** 3 Soi Udomauk 41, Sukhumvit Road,  
Bangchack, Prakhnong, Bangkok 10260

Page 1 of 5

**Equipment:** pH Meter  
**Manufacturer:** METTLER TOLEDO  
**Model:** SevenEasy pH  
**Serial No.:** 1230525212  
**ID No.:** UAE.WAS.003/2553  
**Order No.:** 2302181  
**Operation No.:** 2302181-001  
**Date of Receipt:** 14 March 2023  
**Date of Calibration:** 24 March 2023

**Calibrated by** Mr.Pheraphat Tuanjit **Approved by** [Signature]  
Scientist ( Mr.Nuttapol Niyomchart )  
Specialist, Division of Calibration Laboratory  
Responsible for the Technical Management Team

**Date of Issue:** 24 March 2023

The uncertainties are for a confidence probability of approximately 95%.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme, which has assessed the measurement capability of the laboratory and its feasibility to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65

เอกสารไม่ควบคุม



## Calibration Report

**Certificate No.:** 2302181-001-01  
**Equipment:** pH Meter  
Resolution: 0.01 pH ; 1 mV  
Manufacturer: METTLER TOLEDO Model: SevenEasy pH  
Serial No.: 1230525212 Type: Bench top  
ID No.: UAE.WAS.003/2553

**Date of Calibration:** 24 March 2023 Page 2 of 5

**Location:** Chemical Calibration Laboratory, National Food Institute  
**Environment Condition:** Ambient Temperature: ( 23.4 ± 1.5 ) °C Relative Humidity: ( 52 ± 3 ) %  
**Condition of Equipment:** Good Condition

### Condition of this Results of Calibration

1. Calibration Method In house method : W-CC-002 based on direct measurement by using standard voltage calibrator and certified reference material (CRM)

2. Reference Standards / Certified Reference Material

Instruments	Serial / ID No.	Manufacturer	Certificate No.	Due Date
2.1 DC Voltage Calibrator	2709007	Fuke	22E1959	17 June 2023
2.2 Digital Thermometer	2709007	Fuke	OC-450557-01	30 October 2023
2.3 Thermo-Hygro Meter	NFLBTH053/17	PONPE	TE 650555-01	21 September 2023
Certified Reference Material	Lot No.	Manufacturer	Ref No	Expiry Date
2.4 pH buffer 4.008 (Primary pH buffer Solution)	873608	CPAchem	PH4216.L5	16 February 2025
2.5 pH buffer 6.865 (Primary pH buffer Solution)	873609	CPAchem	PH4217.L5	16 February 2025
2.6 pH buffer 10.01 (Primary pH buffer Solution)	873611	CPAchem	PH4220.L5	16 February 2024
2.7 pH buffer 7.00 (Standard pH buffer Solution)	873612	CPAchem	PH107.L5	16 February 2024

3. This certification is traceable to The International System of Unit (SI Unit)

- 3.1 Instruments No.2.1 through NSC-TIS-TIS 17025 Laboratory Accreditation of Calibration No.008
- 3.2 Instruments No.2.2 through NSC-TIS-TIS 17025 Laboratory Accreditation of Calibration No.0061
- 3.3 Instruments No.2.3 through NSC-TIS-TIS 17025 Laboratory Accreditation of Calibration No.0061
- 3.4 Certified Reference Material No. 2.4 to 2.6 traceable to Primary measurement method- Harned cell using calibrated thermometer, barometer, and nanovoltmeter. The Standard Solution preparation and certified by CPAchem Ltd is accredited to ISO 17034 and ISO/IEC 17025
- 3.5 Certified Reference Material No.2.7 traceable to BSM ReHi H-13 LoPH 25.05.2022; BSM ReHi H-16 LoPH 02.06.2022; BSM ReHi H-13 LoPH 25.05.2022; BSM ReHi H-16 LoPH 02.06.2022, the Standard Solution preparation and certified by CPAchem Ltd is accredited to ISO 17034 and ISO/IEC 17025

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

F-CS-012 Revision: 01 Date: 20-04-65

เอกสารไม่ควบคุม



## Calibration Report

**Certificate No.:** 2302181-001-01  
**Equipment:** pH Meter  
Resolution: 0.01 pH ; 1 mV  
Manufacturer: METTLER TOLEDO Model: SevenEasy pH  
Serial No.: 1230525212 Type: Bench top  
ID No.: UAE.WAS.003/2553

**Date of Calibration:** 24 March 2023 Page 3 of 5

**Calibration Results:**  
1. Calibration of pH Meter ( Manual Temperature Compensation at 25 °C )

Nominal pH	DC Voltage Standard ( mV )	Average Indicator Reading		Uncertainty ( ± mV )	Coverage Factor ( k )
		mV	pH		
0	414.120	414	0.00	0.58	2.00
2	296.614	296	2.00	0.58	2.00
4	177.464	178	4.00	0.58	2.00
6	59.160	59	6.00	0.58	2.00
7	0.000	0	7.00	0.58	2.00
8	-59.158	-59	8.00	0.58	2.00
10	-177.460	-177	10.00	0.58	2.00
12	-296.611	-296	12.00	0.58	2.00
14	-414.117	-414	14.00	0.58	2.00

2. Calibration of pH Meter with Electrode ( Manual Temperature Compensation at 25 °C )

**Equipment:** pH Electrode **Type:** Combined Electrode  
**Manufacturer:** METTLER TOLEDO **Model:** InLab Solids  
**Serial No.:** 1156853 **ID No.:** N/A

**Performance of Electrode system** (Three-Point Calibration at pH 4, pH 7 and pH 10)

Certified Value (25 °C (pH))	Average Indicator Reading		Relative Slope (%)	Uncertainty ( ± pH )	Coverage Factor ( k )
	pH	mV			
4.008	4.01	187	-	0.0075	2.00
6.865	6.86	22	97.86	0.0075	2.00
10.010	10.01	-160	97.66	0.0066	2.00
9.985	9.99	14	-	0.0063	2.00

F-CS-012 Revision: 01 Date: 20-04-65

เอกสารไม่ควบคุม










Equipment : Electronic Balance  
Condition As-Received : Used Item  
Reference : 2304-0459OC-2

Cert.No.: 23MM113  
Page: 3 of 3

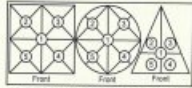

Maximum difference between  
off-center and central loading  
(g)  
0.0001

## 2. Effect of off center loading

A mass of 100 g was placed to various position on the pan.  
The weighing machine reading error obtained is given in the table

Position 1 (g)	Position 2 (g)	Position 3 (g)	Position 4 (g)	Position 5 (g)
-0.0001	-0.0001	0.0000	-0.0001	-0.0001

## 3. Departure from nominal value

Applied Weight (g)	Balance Reading (g)	Correction (g)	Measurement Uncertainty ( $\pm$ mg)	Coverage Factor (k)
Unload	0.00000	0.00000	0.014	2.11
0.05	0.04999	+0.00001	0.015	2.09
0.1	0.09999	+0.00001	0.015	2.07
1	1.00000	0.00000	0.018	2.04
5	5.00000	0.00000	0.026	2.00
20	20.00002	-0.00002	0.045	2.00
50	50.00002	-0.00002	0.080	2.00
80	80.00002	-0.00002	0.15	2.00
100	100.00000	0.00000	0.17	2.00
150	150.00000	0.00000	0.29	2.00
200	199.99999	+0.00001	0.29	2.00

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor  $k$ , providing a level of confidence of approximately 95 %.

-000-

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES & EQUIPMENT CALIBRATION AND TESTING SERVICES  
504/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG RANGKOK 10250  
TEL: 0-2717-3000-27 FAX: 0-2719-9484



Cert. No.: 22TM1490  
Page: 1 of 3

## Certificate of Calibration

Equipment : Hot Air Oven  
Manufacturer : Memmert  
Model : UF 55  
Serial No. : B216.1666  
ID No. : UAE.WAO.027/2559  
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.  
3 Sol Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260  
Location : Lab Floor 2  
Received Order : 19 October 2022  
Calibration Date : 19 October 2022  
Ambient Temperature : ( 26  $\pm$  10 ) °C  
Relative Humidity : ( 50  $\pm$  30 ) %

Calibrated by : Preecha Hishib

Approved by :

Approved Signatory

( ) Pornthippa Tameyakul  
( ) Malee Butkusa  
(✓) Suwit Imjai

Issue Date : 31 October 2022

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services & Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม

A 0046800



Equipment : Hot Air Oven  
Condition As-Received : Used Item  
Reference : 2210-0575OC-1

Cert. No.: 22TM1490  
Page: 2 of 3

## Procedure Used :-

Calibration were conducted using calibration procedure CP-QT02 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD) and Thermocouple Type T.

The temperature scale used was based on ITS-90.

## Condition of this result of calibration

### 1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1) Data Acquisition	34970A	MY41021843	22LM4	10 Jan 2023


2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

Result of Calibration :- ( \* ) Without Adjustment

Function of UUC\* : Temperature Source

Fresh air setting : Close

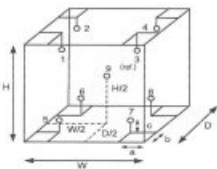

Probe Installation Details :  
a = 5.0 cm  
b = 5.0 cm  
c = 5.0 cm  
Dimension of Chamber :  
D = 0.33 m  
W = 0.40 m  
H = 0.40 m  
Capacity = 0.053 m<sup>3</sup>

Environment during calibration		
	Beginning	Finished
Temp. ( °C )	29	30
REL.Humid. ( % )	47	40
AC Supply ( Volt )	221	220

Ref. Std. ID No.: @ Calibration Point		
Position :	( 104 ) °C	( 140,180 ) °C
1	18-04RTD-01	21-04TC-01
2	18-04RTD-02	21-04TC-02
3	18-04RTD-03	21-04TC-03
4	18-04RTD-04	21-04TC-04
5	18-04RTD-05	21-04TC-05
6	18-04RTD-06	21-04TC-06
7	18-04RTD-07	21-04TC-07
8	18-04RTD-08	21-04TC-08
9 (ref.)	18-04RTD-09	21-04TC-09

เอกสารไม่ควบคุม

a 1133252



Equipment : Hot Air Oven  
Condition As-Received : Used Item  
Reference : 2210-0575OC-1  
Result of Calibration :- ( \* ) Without Adjustment

Function of UUC\* : Temperature Source

Fresh air setting : Close

Calibration Point ( °C )	UUC* Setting ( °C )	UUC* Reading ( °C )	Temperature stability ( $\pm$ °C )	Temperature uniformity ( °C )	Overall Variation ( °C )	Uncertainty ( $\pm$ °C )	Coverage Factor k
104.0	104.0	104.0	0.061	1.3	1.7	0.42	2
140.0	140.0	140.0	0.14	2.3	2.4	1.1	2
180.0	180.0	180.0	0.21	3.5	3.6	1.3	2

Calibration Point ( °C )	Measured Temperature ( °C )								
	1	2	3	4	5	6	7	8	9 (ref.)
104.0	103.076	103.876	103.777	104.124	104.667	104.426	104.012	103.928	104.370
140.0	138.199	139.189	138.808	139.550	140.266	139.622	139.293	139.365	140.369
180.0	177.930	179.267	178.643	179.753	181.011	180.093	179.496	179.743	181.278

Average\* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation

UUC\* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor  $k$ , providing a level of confidence of approximately 95 %.

-000-

เอกสารไม่ควบคุม

a 1133251



## Calibration Certificate

Certificate No.: 2302827-001-01  
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO., LTD.  
Address: 3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong, Bangkok 10260

Page 1 of 4

Equipment: Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR204

Serial No.: C117635043

ID No.: UAE.WAS.012/2564

Order No.: 2302827

Operation No.: 2302827-001

Date of Receipt: 10 May 2023

Date of Calibration: 10 May 2023

Calibrated by Mr. Manas Somsak Specialist  
Approved by (Manasrattana Tungsit)  
Manager, Division of Calibration Laboratory  
Responsible for the Technical Management Team

The uncertainties are for a confidence probability of approximately 95%

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65

## Calibration Report

Certificate No.: 2302827-001-01  
Equipment: Electronic Balance  
Model: XSR204  
Serial No.: C117635043  
Capacity: 220 g  
Manufacturer: METTLER TOLEDO  
Resolution: 0.0001 g  
ID No.: UAE.WAS.012/2564

Page 2 of 4

Date of Calibration: 10 May 2023  
Environment Condition: Ambient Temperature: 23.4 ± 0.2 °C Relative Humidity: 43.4 ± 0.9 %

Place of Calibration: Balance room (Water Analysis Unit), UNITED ANALYST AND ENGINEERING CONSULTANT CO., LTD.

Condition of Equipment: Good Condition

Condition of This Results of Calibration:

1. Calibration Method: NFI Method W-MA-001 In-House Method based on UKAS Lab 14: 2019

2. Reference Standards:

Reference Standard	Model	Serial No.	Calibrated By	Certificate No.	Due Date
Standard Weight Class E2	1mg to 200g	B50567572	TCS	M23040535	8 April 2024

Instrument	Model	Serial No.	Calibrated By	Certificate No.	Due Date
Thermo-Hygro Meter	608-H1	NFI.BTH 016/23	Quality Reborn	QR23-0489	21 February 2024

3. This certification is traceable to SI UNIT

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

Calibration Results:


1. Repeatability of Reading:

Nominal Value (g)	Standard Deviation of Reading (g)
100	0.000032
200	0.000032

2. Off-Center Error:

A mass of 100 g was placed and moved to various position on pan.

The balance reading obtained is given in the table.

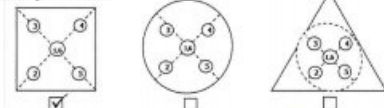

1	2	3	4	5	6	(Maximum Difference)
(g)	(g)	(g)	(g)	(g)	(g)	(g)
100.0002	100.0002	100.0002	100.0002	100.0003	100.0002	0.0001

F-CS-012 Revision: 01 Date: 20-04-65

## Calibration Report

Certificate No.: 2302827-001-01  
Equipment: Electronic Balance  
Model: XSR204  
Serial No.: C117635043  
Capacity: 220 g  
Manufacturer: METTLER TOLEDO  
Resolution: 0.0001 g  
ID No.: UAE.WAS.012/2564

Page 3 of 4

Date of Calibration: 10 May 2023

Calibration Results: (Continued)

Calibration Range: 0 - 200 g

Calibration Adjustment: Internal Calibration

3. Departure from Nominal Value:

Nominal Value (g)	Standard Value (g)	Average Reading (g)	Correction (g)	Uncertainty (± g)	Coverage Factor k
Unread	0.00000	0.0000	0.0000	0.000005	2.00
0.01	0.01000	0.0100	0.0000	0.000005	2.00
0.02	0.02001	0.0200	0.0000	0.000005	2.00
0.05	0.05000	0.0500	0.0000	0.000005	2.00
0.1	0.10001	0.1000	0.0000	0.000005	2.00
0.2	0.20001	0.2000	0.0000	0.000005	2.00
0.5	0.50002	0.5000	0.0000	0.000005	2.00
1	1.00000	1.0000	0.0000	0.000005	2.00
2	2.00002	2.0000	0.0000	0.000005	2.00
3	3.00003	3.0000	0.0000	0.000007	2.00
5	5.00002	5.0000	0.0000	0.000007	2.00
10	10.00003	10.0000	0.0000	0.000008	2.00
20	20.00003	20.0000	0.0000	0.000009	2.00
30	30.00004	30.0000	0.0000	0.000009	2.00
40	40.00007	40.0000	0.0000	0.000011	2.00
45	45.00009	45.0001	0.0000	0.000013	2.00

F-CS-012 Revision: 01 Date: 20-04-65

## Calibration Report

Certificate No.: 2302827-001-01  
Equipment: Electronic Balance  
Model: XSR204  
Serial No.: C117635043  
Capacity: 220 g  
Manufacturer: METTLER TOLEDO  
Resolution: 0.0001 g  
ID No.: UAE.WAS.012/2564

Page 4 of 4

Date of Calibration: 10 May 2023

Calibration Results: (Continued)

Calibration Range: 0 - 200 g

Calibration Adjustment: Internal Calibration

3. Departure from Nominal Value:

Nominal Value (g)	Standard Value (g)	Average Reading (g)	Correction (g)	Uncertainty (± g)	Coverage Factor k
50	50.00003	50.0000	0.0000	0.000011	2.00
55	55.00005	55.0000	0.0000	0.000012	2.00
60	60.00004	60.0000	0.0000	0.000012	2.00
65	65.00005	65.0000	0.0000	0.000013	2.00
70	70.00006	70.0001	-0.0001	0.000013	2.00
75	75.00008	75.0002	-0.0001	0.000013	2.00
80	80.00007	80.0002	-0.0001	0.000014	2.00
85	85.00009	85.0002	-0.0001	0.000014	2.00
90	90.00010	90.0002	-0.0001	0.000015	2.00
100	100.00006	100.0002	-0.0001	0.000016	2.00
120	120.00009	120.0002	-0.0001	0.000018	2.00
150	150.00009	150.0002	-0.0001	0.000021	2.00
200	200.00016	200.0003	-0.0001	0.000028	2.00

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k providing a level of confidence of approximately 95 %.

F-CS-012 Revision: 01 Date: 20-04-65





Cert. No.: 23TM249  
Page : 1 of 3

## Certificate of Calibration

**Equipment :** BOD Incubator  
**Manufacturer :** Arco  
**Model :** UC4-1320  
**Serial No. :** 13URC4S013201  
**ID No. :** UAE.WAO.015/2561  
**Submitted by :** United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260  
**Location :** Lab Floor 2  
**Received Order :** 15 February 2023  
**Calibration Date :** 15 February 2023  
**Ambient Temperature :** ( 26 ± 10 ) °C  
**Relative Humidity :** ( 50 ± 30 ) %  
**Calibrated by :** Preecha Hiahib

**Approved by :**  
( ) Pornthippa Tamayakul  
(✓) Malee Butkruea  
( ) Suwit Imjai

**Issue Date :** 24 February 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3: Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม

A 0051476



**Equipment :** BOD Incubator  
**Condition As-Received :** Used Item  
**Reference :** 2302-0297OC-1  
**Procedure Used :-**

Cert. No.: 23TM249  
Page : 2 of 3

Calibration were conducted using calibration procedure CP-OT02 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector ( RTD ).  
The temperature scale used was based on ITS-90.

### Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1 ) Data Acquisition	34972A	MY57013711	22LM93	02 Jul 2023


2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

**Result of Calibration :-** ( \* ) Without Adjustment

**Function of UUC\* :** Temperature Source

**Fresh air setting :** Not Available

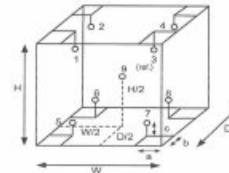

### Probe Installation Details :

a = 10 cm  
b = 10 cm  
c = 10 cm

### Dimension of Chamber :

D = 0.62 m  
W = 1.2 m  
H = 1.2 m  
Capacity = 0.89 m<sup>3</sup>

Environment during calibration		
	Beginning	Finished
Temp. ( °C )	29	31
REL.Humid. ( % )	63	67
AC Supply ( Volt )	220	220

Position :	Ref. Std. ID No.:
1	22-18RTD-2/1
2	18RTD-2/2
3	18RTD-2/3
4	18RTD-2/4
5	18RTD-2/5
6	18RTD-2/6
7	18RTD-2/7
8	18RTD-2/8
9 (ref.)	18RTD-2/9

เอกสารไม่ควบคุม

0 4440547



**Equipment :** BOD Incubator  
**Condition As-Received :** Used Item  
**Reference :** 2302-0297OC-1  
**Result of Calibration :-** ( \* ) Without Adjustment  
**Function of UUC\* :** Temperature Source  
**Fresh air setting :** Not Available

Cert. No.: 23TM249  
Page : 3 of 3

Calibration Point ( °C )	UUC* Setting ( °C )	UUC* Reading ( °C )	Temperature stability ( ± °C )	Temperature uniformity ( °C )	Overall Variation ( °C )	Uncertainty ( ± °C )	Coverage Factor k
20.0	20.0	19.3	0.32	0.57	1.0	0.60	2

Calibration Point ( °C )	Measured Temperature ( °C )								
	Position								
20.0	1	2	3	4	5	6	7	8	9 (ref.)
	20.086	19.916	20.386	19.976	19.973	19.836	19.837	19.821	19.949

**Average\* :** The average of 30 values in each position.

**Temperature stability :** One-half of the greatest maximum difference of measured temperature at any one sensor.

**Temperature uniformity :** The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

**Overall Variation :** The Difference of the maximum and minimum measured temperatures throughout observation.

**UUC\* :** Unit Under Calibration

**Note :** The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-000-

เอกสารไม่ควบคุม

a 1149512



Cert. No.: 23TM375  
Page : 1 of 3

## Certificate of Calibration

**Equipment :** BOD Incubator

**Manufacturer :** ARCO

**Model :** UR-1320

**Serial No. :** -

**ID No. :** UAE.WAO.018/2551

**Submitted by :** United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260

**Location :** Lab Floor 2

**Received Order :** 11 April 2023

**Calibration Date :** 12 April 2023

**Ambient Temperature :** ( 26 ± 10 ) °C

**Relative Humidity :** ( 50 ± 30 ) %

**Calibrated by :** Krisda Malee

**Approved by :**  
( ) Pornthippa Tamayakul  
(✓) Malee Butkruea  
( ) Suwit Imjai

**Issue Date :** 24 April 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3: Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม

A 0053360



Equipment : BOD Incubator  
Condition As-Received : Used Item  
Reference : 2304-0156OC-2  
Procedure Used :-

Cert. No.: 23TM375  
Page : 2 of 3

Calibration were conducted using calibration procedure CP-OT02 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector ( RTD ).

The temperature scale used was based on ITS-90.

#### Condition of this result of calibration

##### 1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1 ) Data Acquisition	34972A	MY59003411	22LM165	26 Nov 2023


2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

Result of Calibration :- ( \* ) Without Adjustment

Function of UUC\* : Temperature Source

Fresh air setting : Not Available

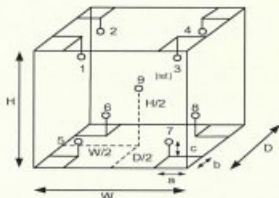

#### Probe Installation Details :

a = 10 cm  
b = 10 cm  
c = 10 cm

#### Dimension of Chamber :

D = 0.62 m  
W = 1.2 m  
H = 1.2 m  
Capacity = 0.89 m<sup>3</sup>

Environment during calibration		
	Beginning	Finished
Temp. ( °C )	26	27
REL.Humid. ( % )	42	45
AC Supply ( Volt )	219	220

Position :	Ref. Std. ID No.:
1	20RTD-2/1
2	20RTD-2/2
3	20RTD-2/3
4	20RTD-2/4
5	20RTD-2/5
6	20RTD-2/6
7	20RTD-2/7
8	20RTD-2/8
9 (ref.)	20RTD-2/9

เอกสารไม่ควบคุม

a 1158259



Equipment : BOD Incubator  
Condition As-Received : Used Item  
Reference : 2304-0156OC-2

Cert. No.: 23TM375  
Page : 3 of 3

Result of Calibration :- ( \* ) Without Adjustment

Function of UUC\* : Temperature Source

Fresh air setting : Not Available

Calibration Point ( °C )	UUC* Setting ( °C )	UUC* Reading ( °C )	Temperature stability ( ± °C )	Temperature uniformity ( °C )	Overall Variation ( °C )	Coverage Factor k
20.0	20.0	20.0	0.48	0.42	1.2	2

Calibration Point ( °C )	Measured Temperature ( °C )									Uncertainty ( ± °C )
	1	2	3	4	5	6	7	8	9 (ref.)	
20.0	20.040	20.170	20.263	20.093	19.749	19.704	19.920	20.191	20.020	0.66

Average\* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC\* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่ควบคุม

a 1158259



จุฬาลงกรณ์มหาวิทยาลัย  
ศูนย์บริการห้องปฏิบัติการ  
Foundation for Industrial Development National Food Institute  
Food Industrial Laboratory Service Center

## Verification Certificate

Certificate No.: 2302413-001-01  
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO., LTD.  
Address: 3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchack, Prakanong, Bangkok 10260

Page 1 of 4

Equipment: HEATING BLOCK DIGESTION  
Manufacturer: FOSS  
Model: 2520  
Serial No.: 91794469  
ID No.: UAE.WAS.011/2560  
Order No.: 2302413  
Operation No.: 2302413-001  
Date of Receipt: 28 March 2023  
Date of Calibration: 30-31 March 2023

Calibrated by Mr.Nuttapol Niyomchat Specialist  
Approved by ( Mr.Pheraphat Tuanjit )  
Manager, Division of Calibration Laboratory  
Responsible for the Technical Management Team  
Date of Issue: 10 April 2023

The uncertainties are for a confidence probability of approximately 95 %.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65

เอกสารไม่ควบคุม



จุฬาลงกรณ์มหาวิทยาลัย  
ศูนย์บริการห้องปฏิบัติการ  
Foundation for Industrial Development National Food Institute  
Food Industrial Laboratory Service Center

## Verification Report

Certificate No.: 2302413-001-01  
Equipment: HEATING BLOCK DIGESTION  
Model: 2520 Serial No.: 91794469  
Resolution: 1 °C ID No.: UAE.WAS.011/2560  
Manufacturer: FOSS

Date of Calibration: 30-31 March 2023 Page 2 of 4

Location: Laboratory Room, NATIONAL FOOD INSTITUTE  
Environment Condition: Ambient Temperature ( 25 ± 3 ) °C  
Relative Humidity ( 55 ± 15 ) %  
Line Voltage ( 220 ± 10 ) Volt

#### Condition of this results of Calibration:

1. This instrument was calibrated by insert standard thermocouples type R into its heating block digestion and compared to temperature obtained from reference standards thermometer at calibrated point.

- The temperature scale used was based on ITS - 90 .

- All data show below were final values and the initial data may be obtained upon request.

#### 2. Reference Standard Instrument :

Instrument	Model	Serial No.	Certificate No.	Due Date	Through
Digital Thermometer with Thermocouple	34970A	MY4465536 / MY45304652	TC22/0044	5-May-2023	N.M. Technical Center Laboratory
	Type R	TC910A-102 / CH9101-102			

3. This certificate is traceable to international system of units (SI Units).

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

6. Condition of Calibrated item : Good

UUC\* Description

Time of Record - Hour 30 Minute At 380 °C

7. Result of Calibration : ☒ Without adjustment ☐ After adjustment

F-CS-009 Revision: 01 Date: 20-04-65

เอกสารไม่ควบคุม





## Verification Report

Certificate No.: 2302413-001-01  
Equipment: HEATING BLOCK DIGESTION  
Model: 2520 Serial No.: 91794469  
Resolution: 1 °C ID No.: UAE.WAS.011/2560  
Manufacturer: FOSS  
Date of Calibration: 30-31 March 2023  
Calibration point: 380 °C  
Calibration result: Page 3 of 4

Block No.	UUC* Setting (°C)	UUC* Reading (°C)	Stability (±°C)	Standard Thermometer (°C)	Uncertainty (±°C)
1	380	380	0.96	377.74	2.1
2	380	380	0.40	377.28	2.1
3	380	380	1.18	377.82	2.1
4	380	380	0.44	377.19	1.6
5	380	380	0.11	377.30	1.6
6	380	380	0.14	377.90	1.6
7	380	380	1.17	373.85	2.1
8	380	380	0.33	376.96	2.1
9	380	380	0.14	374.18	2.1
10	380	380	0.96	378.56	2.0
11	380	380	1.04	378.34	2.0
12	380	380	0.35	378.06	2.0
13	380	380	0.48	377.05	1.6
14	380	380	0.38	379.19	1.6
15	380	380	0.50	377.48	1.6
16	380	380	0.48	378.33	1.7
17	380	380	0.71	377.60	1.7
18	380	380	0.35	376.77	1.7
19	380	380	0.84	377.06	1.8
20	380	380	0.41	378.58	1.8


Note:  
- UUC\* = Unit Under Calibration  
- Immersion depth of standard thermometer in tube level high of sand is equal heater plate of UUC.  
- Stability = One-half of the greatest maximum difference of measured temperatures at one sensors, for at least half an hour after reaching steady state.

FCS-009 Revision: 01 Date: 20-04-65

## Verification Report

Certificate No.: 2302413-001-01  
Equipment: HEATING BLOCK DIGESTION  
Model: 2520 Serial No.: 91794469  
Resolution: 1 °C ID No.: UAE.WAS.011/2560  
Manufacturer: FOSS  
Date of Calibration: 30-31 March 2023  
Calibration point: 380 °C  
Calibration result: Continued Page 4 of 4

Figure 1. Location of Reference Standard and Block Diagram of Digestion Unit

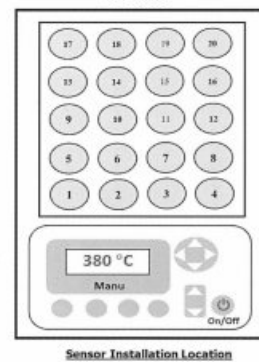

Note:  
- UUC\* = Unit Under Calibration  
- Immersion depth of standard thermometer in tube level high of sand is equal heater plate of UUC.  
- Stability = One-half of the greatest maximum difference of measured temperatures at one sensors, for at least half an hour after reaching steady state.

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor  $k=2$  providing a level of confidence of approximately 95 %.

\*\*\*\*\* End \*\*\*\*\*

FCS-009 Revision: 01 Date: 20-04-65

**FOSS**  
Customer Service Report

FOSS South East Asia  
3388 Srinrat Building, 25th - 26th Floor, Unit No. 3388/90,  
Rama IV Road, Klongton - Klongtoey, Bangkok, Thailand 10110

Report No.: 6623

Date: 30/03/23  
Customer: United Analyst and Engineering  
Address: 115/115  
Instrument: KT5100  
Serial: 318890052

Hours: Travel To Customer: 3.00 - 3.30  
Start: 3.00  
Finish: 3.30  
Labour: 3.00 - 4.00  
Travel From Customer: 4.00 - 4.30  
Start: 4.00  
Finish: 4.30

Application	Special	Job Type	Standard	Training
Normal	Courtesy Visit	Installation	X	
Distributor	PMA Onboarding	Quote		
Internal	Warranty	Repair		
Digital Service	Sales Support	Remote		

PO/Quote Number: 6623

PMA Type: PMA Type Contract No. Contract No.

Details of Work / Test

Details of Work / Test	Condition / Status
Unpack, inspect and test the instrument	OK
Connect the instrument to the computer	OK
Install the software	OK
Calibrate the instrument	OK
Test the instrument	OK
Return the instrument	OK

Instrument Ready for Use: ☒ OK ☐ Not OK

Part No.	Batch	Description	Qty

I confirm this report is accurate and complete

Signed FOSS: [Signature] Signed Customer: [Signature]  
Name: [Name] Name: [Name]

Would you be willing to participate in a brief survey in order to tell us how we performed? ☐

เอกสารไม่ควบคุม

**FOSS**  
Customer Service Report

FOSS South East Asia  
3388 Srinrat Building, 25th - 26th Floor, Unit No. 3388/90,  
Rama IV Road, Klongton - Klongtoey, Bangkok, Thailand 10110

Report No.: 6534

Date: 28/12/22  
Customer: United Analyst and Engineering  
Address: 115/115  
Instrument: KT5100  
Serial: 318890052

Hours: Travel To Customer: 3.00 - 3.30  
Start: 3.00  
Finish: 3.30  
Labour: 3.00 - 4.00  
Travel From Customer: 4.00 - 4.30  
Start: 4.00  
Finish: 4.30

Application	Special	Job Type	Standard	Training
Normal	Courtesy Visit	Installation	X	
Distributor	PMA Onboarding	Quote		
Internal	Warranty	Repair		
Digital Service	Sales Support	Remote		

PO/Quote Number: 6534

PMA Type: PMA Type Contract No. Contract No.

Details of Work / Test

Details of Work / Test	Condition / Status
Unpack, inspect and test the instrument	OK
Connect the instrument to the computer	OK
Install the software	OK
Calibrate the instrument	OK
Test the instrument	OK
Return the instrument	OK

Instrument Ready for Use: ☒ OK ☐ Not OK

Part No.	Batch	Description	Qty

I confirm this report is accurate and complete

Signed FOSS: [Signature] Signed Customer: [Signature]  
Name: [Name] Name: [Name]

Would you be willing to participate in a brief survey in order to tell us how we performed? ☐

เอกสารไม่ควบคุม



## Kjeltec™ 8100 Distillation Unit

This IQ applies to Kjeltec™ 8100 Distillation Unit manufactured by FOSS Analytical. The installation is performed by FOSS trained service personnel.

## 1 Intended Use

Kjeltec 8100 is intended for laboratory use analyzing parameters as specified in FOSS Analytical AB's Application Notes.

## 2 Purpose

This installation Qualification is designed to assure that:

- The Kjeltec instrument is received complete, with all required parts in good condition.
- The location of the instrument is environmentally and ergonomically suitable
- The instrument is assembled and configured correctly
- Suitable electricity and water are supplied to the instrument, see table 2 for requirements.

## 3 Identification

Description	Serial Number
Kjeltec 8100 Distillation Unit	918 89052

## Dedicated Analytical Solutions

FOSS Analytical AS  
69 Slangerupgade  
DK-3400 Hillerød  
Denmark

Tel +45 7010 3370  
Fax +45 7010 3371  
E-mail support@foss.dk  
Web www.foss.dk

FOSS Analytical AB  
Box 70  
SE-263 21 Högabö  
Sweden

Tel +46 42 361500  
Fax +46 42 340349  
E-mail support@foss.se  
Web www.foss.se

Customer Support, 6903 7242 / Rev. 1

เอกสารไม่ควบคุม

190

## 4 Control of Received Equipment

## 4.1 Verify that the correct instrument type and accessory kit items are received and in proper condition

The packing list (shipped with the instrument) specifies all the items. The installer will verify that all items are received as shipped on the packing list. For each item listed, verify that the acceptance criteria are met. If so, write "Y" in the right column of the table immediately following.

Packing List Item	Acceptance Criteria	Pass (Y/N)
Kjeltec 8100 Distillation Unit	No visible damage, received in undamaged FOSS Analytical's standard shipping container	Y
Accessory kit, according to packing list	Included. No visible damage, received in undamaged FOSS Analytical's standard shipping container	Y
Handling device for digestion tube	Included. No visible damage.	Y
Tanks with level sensors for Waste, Alkali and Water	Included. No visible damage.	Y
Receiver flask	Included. No visible damage.	Y
One digestion tube 250ml	Included. No visible damage.	Y
One digestion tube 100 ml	Included. No visible damage.	Y
Tube adapter	Included. No visible damage.	Y
User manual	Kjeltec 8100 Distillation Unit	Y
Owners guide	Kjeltec 8100 Distillation Unit	Y
Quick guide	Kjeltec 8100 Distillation Unit	Y
Spare parts manual	Kjeltec 8100 Distillation Unit	Y
Application notes	AN 300 included AN 303 included	Y

Customer Support, 6903 7242 / Rev. 1

เอกสารไม่ควบคุม

200

## 5 Installation

## 5.1 The equipment must be installed in a suitable location with power, water and draining available

Verify that the instrument installation site meets the acceptance criteria given in the table below. If so, write "Y" in the right column of the table immediately following.

Location Requirements	Acceptance Criteria	Pass (Y/N)
Adequate space for instrument	Dimensions 48x58x69 cm	Y
AC supply available for instrument	200-240 V 50/60Hz	Y
Current	10 A	Y
Cold water supply available	2 l/min at 30°C	Y
Drain	For cooling water and waste (depending on local waste disposal legislation)	Y
Ambient temperature	Max. 40°C	Y
Ambient humidity	Max. 80% relative	Y
Internal fuses	T10A AH	Y


Customer Support, 6903 7242 / Rev. 1

เอกสารไม่ควบคุม

200

## 5.2 The instrument must be assembled correctly

Verify that all tubes are correct connected. If so, write "Y" in the right column of the table immediately following.

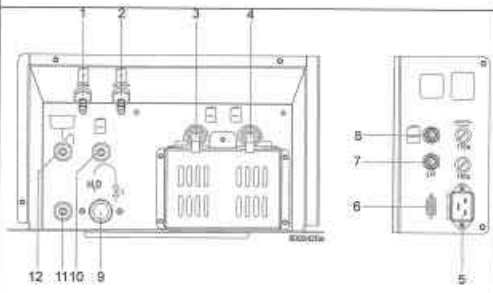Instrument Tubing Connections	Acceptance Criteria	Pass (Y/N)
 <p>1. Deionised water in (steam generator) 2. Deionised water in (dilution water) 3. *) Receiver solution in 4. Alkali in 5. Power 6. Not used 7. External titration module 8. Level sensors 9. Cooling water in (tap water) 10. Waste water out (tube drain vessel) 11. Drain 12. Cooling water out (tap water) *) Only on Kjeltec 8200</p>	Visual verification by installer	Y

Customer Support, 6903 7242 / Rev. 1

เอกสารไม่ควบคุม

400

### 5.3 The instrument should be assembled and powered up

Connect the distilling unit to the power supply. Perform the start up procedure and check that the expected response is obtained. If so, write "Y" in the right column of the table immediately following.

Action	Expected Response	Pass (Y/N)
Switch on the power	The instruments start up and the self test will run. The sample counter shows the number of analysed samples since first power and the Software Version shows the version of the instruments software.	Y
	After start-up, Program 1 is loaded and the Analyse menu is displayed.	Y
Turn on the cold water tap	No visible reaction	Y
Press the "Manual" view	The Manual menu is opened	Y
Open the door with the handle, place the test tube and receiver flask in position. Close the door.		Y
Select <b>Dilution</b> and press <b>Start</b>	Water is added to the tube	Y
Select <b>Alkali</b> and press <b>Start</b>	Alkali is added to the tube	Y
Select <b>Steam</b> and press start	After heating up, steam is entering the tube	Y
Select <b>Drain</b> and press <b>Start</b>	The tube is drained	Y

## 6 Summary of Deviations/Comments

Deviations from above requirements are specified below and any corrective actions are noted.

Deviation	Action	Comment

## 7 IQ Documentation

Upon successful completion and recording of all instructions above, sign and date this sheet below. If required by customer, leave one signed copy with instrument.

If customer's internal procedures require further reporting or witnessing of results, execute those procedures as required.

Installed By: Pannipa Onnon  
 Company: Foss SEA  
 Customer Name: United Analyst and Engineering  
 Company: United Analyst and Engineering  
 Date completed: July 25, 2022

FOSS

Operation Qualification

## Kjeltec™ 8100 Distillation Unit

This OQ applies to Kjeltec 8100 Distillation Unit manufactured by FOSS Analytical. The operation qualification is performed by FOSS trained service personnel.

### 1 Intended Use

Kjeltec 8100 is intended for laboratory use analyzing parameters as specified in FOSS Analytical Application Notes.

### 2 Purpose

This procedure is designed to test the function of the instrument according to factory test specifications:

- Alkali volume
- Distillation Accuracy
- Distillation Repeatability

### 3 Identification

Description	Serial Number
Kjeltec 8100 Distillation Unit, 200-240 V 50/60 Hz	91889052

## 4 Performance

### 4.1 Verify the dispensed volumes of reagents

Note! To verify the dispensed volumes of reagents a triple test should be done to be statistic correct. Then calculate a mean value.

1. Choose "Manual" in the menu. (When starting up the instrument Program 1 is loaded)
2. Open the safety door by pressing **Open** and place a tube in the instrument. Close the safety door.

#### Water

1. Press **Dilution** and then press **Start**. 80 ml of water will be filled into the tube.
2. Measure the collected water in a graduated measuring glass and note the result in table 1 below.
3. Check acceptance criteria in the table and make the judgment if passed or not.

Note! If the water volume needs to be calibrated, go to 4.8.5 Dilution Pump Calibration in the User Manual.

#### Alkali

1. Press **Alkali** and then press **Start**. 50 ml of alkali will be filled into the tube.
2. Measure the collected alkali in a graduated measuring glass and note the result in table 1 below.
3. Check acceptance criteria in the table and make the judgment if passed or not.

Table 1 Volume control

Test	Result	Expected result	Passed (Y/N)
Water volume	73 ml 75 ml 71 ml Mean 73.0 ml	76- 84 ml	Y
Alkali volume	43 ml 42 ml 53 ml Mean 42.7 ml	47- 54 ml	Y

### Dedicated Analytical Solutions

FOSS Analytical A/S  
 69 Slangerupgade  
 DK-3400 Hillerød  
 Denmark  
 Tel +45 7010 3370  
 Fax +45 7010 3371  
 E-mail support@foss.dk  
 Web www.foss.dk

FOSS Analytical AB  
 Box 10  
 SE-283 21 Högabäck  
 Sweden  
 Tel +46 42 361500  
 Fax +46 42 340049  
 E-mail support@foss.dk  
 Web www.foss.dk

## 4.2 Verify the distillation procedure, accuracy and precision

The distillation principle is to convert ammonium ( $\text{NH}_4^+$ ) into ammonia ( $\text{NH}_3$ ) by using an alkali (NaOH) and thereafter steam distil it into a receiver flask containing boric acid and titrate with standard acid solution using colorimetric end-point detection. Ammonium sulphate, a substance with known ammonia content, can be used to check the accuracy of the distillation. The recovery is calculated from obtained result.

The way to perform this test will be described in the following.

### Chemical Check

Use ammonium sulphate ( $\text{NH}_4)_2\text{SO}_4$ , purity > 99.5 % \*)

Mol. weight = 132.14 g/mol, Nitrogen content in ammonium sulphate (99.5 %) = 21.09% \*)

### Analysis conditions according to AN 300

Water	80 ml
Alkali	50 ml NaOH (40%w/w)
Receiver solution	30 ml boric acid (4%)
Distillation time	5 minutes
SAFE	5 seconds
Titrant	0.2N HCl

### For reagent preparation see Appendix A

- Start the instrument and run two blanks without chemicals according to above analysis conditions, distil into a receiver flask containing boric acid. Titrate with a standard acid solution using colorimetric end-point detection. If the blanks are less than 0.2 ml continue with the recovery tests;
- Weigh 0.15 g ammonium sulphate into a tube. Prepare 6 samples (tubes).
- Run the six samples according to above analysis conditions. Titrate with a standard acid solution using colorimetric end-point detection.
- Calculate the recovery according to below equations. Expected results of recovery should be 100%±1%.

Recovery test	Result	Expected result	Passed (Y/N)
Blank value (water blank)	1. <u>0.05</u> ml 2. <u>0.14</u> ml	0.05-0.20 ml	Y
Recovery	1. <u>100.50</u> % 2. <u>100.30</u> % 3. <u>100.45</u> % 4. <u>99.81</u> % 5. <u>99.97</u> % 6. <u>100.81</u> %		
Accuracy	Mean Value: <u>100.0%</u>	99-101 %	Y
Precision	SD: <u>0.557</u>	SD <1 %	Y

Customer Support, 6003 7246 / Rev. 1

เอกสารไม่ควบคุม

307

\*) Note! Please also note that the below calculations must be adjusted if other purity levels of ammonium salts are used. A certificate from the chemical supplier should be available

Purity	Nitrogen content
99,5%	21.09%
99,6%	21.12%
99,7%	21.14%
99,8%	21.16%
99,9%	21.18%

$$\% \text{ Nitrogen} = \frac{(ml_{\text{sample}} - ml_{\text{blank}}) \times N \times 14,007 \times 100}{mg_{\text{sample}}} \quad \begin{matrix} \nearrow 0.1095 \\ 21.72 \end{matrix}$$

N = Normality of titrant to 4 places of decimal.

$$\% \text{ Recovery} = \frac{\% \text{ Nitrogen}}{21.09} \times 100$$

$$mg \text{ sample} = \frac{(ml_{\text{sample}} - ml_{\text{blank}}) \times N \times 14,007 \times 100}{\% \text{ Nitrogen}}$$

① 0.1598 21.56  
②  
③  
④  
⑤  
⑥

Customer Support, 6003 7246 / Rev. 1

เอกสารไม่ควบคุม

417

## 5 Summary of Deviations/Comments

Deviations from above requirements are specified below and any corrective actions are noted.

Deviation	Action	Comment

## 6 OQ Documentation

Upon successful completion of tests above, sign and date this sheet below. If required by customer, leave one signed copy with instrument.

If customer's internal procedures require further reporting or witnessing of results, execute those procedures as required.

Performed By: \_\_\_\_\_  
Company: \_\_\_\_\_  
Customer Name: \_\_\_\_\_  
Company: \_\_\_\_\_  
Date completed: \_\_\_\_\_

## 7 Appendix A

### 7.1 Preparation of Reagents

#### 7.1.1 Alkali

To convert ammonium into ammonia an excess of sodium hydroxide is necessary.

Use 400 g NaOH per litre of solution. Commercially available in concentrations up to 50 %. Do not use concentrations above 40 % as this will lead to crystal formation impairing the function of the pumps. If you can only buy concentrations > 40 %, dilute it before use.

#### 7.1.2 Titrant acid, determination of concentration

To be able to achieve accurate nitrogen / protein results, one must be quite sure that the HCl (hydrochloric acid) concentration is what it is supposed to be. A titration against a predetermined solution of sodium carbonate as described below is thus necessary. Incorrect HCl concentration can otherwise cause substantial errors.

- Standard substance**  
Weigh approx. 10 g of anhydrous sodium carbonate ( $\text{Na}_2\text{CO}_3$ ). Use a mortar to make a fine powder. Dry it for 1 h at 265 °C or 2 h at 200 °C. After cooling in a desiccator, transfer the sodium carbonate to a beaker with a tight lid. Store it in a desiccator.
- Indicator solutions**  
Dissolve 0.1 g methyl red in 100 ml methanol. Dissolve 0.1 g bromocresol green in 100 ml methanol.
- Procedure**  
Weigh approx. 0.4 g of the standard substance, using an analytical balance, note the weight ( $W_1$ ). Transfer the sodium carbonate to a receiver flask and add 40 ml of  $\text{H}_2\text{O}$  (distilled or deionized). Add 8 drops from each of the indicator solutions. Titrate to pink. Note the amount in ml used ( $A_1$ ). Boil this solution for a few minutes. The solution will turn green. Cool rapidly to room temperature under running water. Continue the titration until the next pink colour change occurs. Note also this volume ( $A_2$ ). Boil the solution for a few minutes. Cool rapidly to room temperature under running water. Continue the titration until the next pink colour occurs. Note also this volume ( $A_3$ ).  
**Note!** Temperature changes will influence the volume and the concentration of the titrant solution. The working temperature of the titrant should approximate that of its temperature during standardization. If temperature corrections are necessary, sufficient accuracy may be obtained by use of a correction table. (AOAC 942.25)

Customer Support, 6003 7246 / Rev. 1

เอกสารไม่ควบคุม

307

Customer Support, 6003 7246 / Rev. 1

เอกสารไม่ควบคุม

417



## 7.2 Calculation

$$\text{Molarity (M)} = \frac{18,870 \times W_1}{(A_1 + A_2 + A_3)}$$

**Note!** Concentration must be accurate to four digits, i.e. 0.2000 M.

**Note!** The colour change of this official procedure (AOAC 936.15) may be difficult to see, therefore a pH meter or a mixed indicator (e.g. 0.1 g Methyl red and 0.1 g Bromocresol green in 100 ml methanol) will make it much easier to perform.

## 7.3 Receiver Solution

Boric acid 4 % with bromocresol green / methyl red indicator solution

In order to obtain accurate results the receiver solution is adjusted so that a small (0.05–0.20 ml) positive blank is obtained when running a blank sample. The 4 % boric acid receiver solution is prepared by dissolving 400 g of boric acid in about 5–6 l very hot deionized water. Mix and add more hot deionized water to a volume of about 9 l. Cool the solution to room temperature and add 100 ml of bromocresol green solution (100 mg in 100 ml methanol) and 70 ml of methyl red solution (100 mg in 100 ml of methanol). Dilute to 10 l with deionized water and mix carefully.

**Note!** The addition of alkali is to achieve a positive blank value. This should, however, be kept between 0.05 - 0.20 ml titrant, to obtain good repeatability when testing blanks.

Adjustment of the boric acid is made by the following procedure:

1. Transfer 25 ml boric acid solution to a receiver flask and add 100 ml of distilled water. If the solution in the flask is still red, titrate with 0.1 M sodium hydroxide solution until a neutral grey colour is obtained. Calculate the amount of sodium hydroxide solution necessary to adjust the boric acid solution in the 10 l flask with the formula: ml 1.0 M alkali = ml titrant x 40.
2. Add the calculated amount of 1.0 M alkali solution to the boric acid solution. Mix.
3. To check proceed as follows using 25 ml of the boric acid solution. Run a blank. If the value of this blank is high (0.5 ml of 0.2 M HCl) the boric acid is incorrectly adjusted. This might create irregular blanks. For correction add HCl directly into the boric acid tank, mix it carefully and repeat until a reading of 0.05 - 0.20 ml HCl is obtained. If a positive blank is not achieved, add further small quantities of 1 M NaOH and repeat the check until a satisfactory value is achieved.

## 4.2 Ongoing Qualification Tests

### Block Temperature

The temperature for the digestion is limited by the boiling point for the sulphuric acid, this can be increased by adding a salt ( $K_2SO_4$ ) to the digestion mixture. It's important that the optimal ratio between acid and salt is kept; please follow recommendation in AN 300 or suggested procedures for a specific kind of sample material.

The block temperature itself can be controlled external by inserting a temperature probe in the intended hole in the aluminium block (front row of holes).

Use the reagents and method procedure specified in AN 300. Use only reagents of recognized analytical grade, unless otherwise specified and distilled or demineralised water or water of equivalent purity.

*Suggested standard material for internal quality control:*

Ammonium sulphate  $[(NH_4)_2SO_4]$ , min. 99.5 % (mass fraction), with certified purity.

**Note: The above chemical is usually readily available with a certificate specifying the purity.**

Alternatively ammonium iron(II) sulphate,  $(NH_4)_2 Fe(SO_4)_2 \times 6 H_2O$ , with certified purity may be used.

Tryptophan ( $C_{11}H_{12}N_2O_2$ ), minimum assay 99 % (mass fraction). Nitrogen content 137.2 g/kg. Do not dry in an oven before use.

Acetanilide ( $C_8H_9NO$ ), minimum assay 99 % (mass fraction). Nitrogen content 103.6 g/kg. Do not dry in an oven before use.

Sucrose, ( $C_{12}H_{22}O_{11}$ ), with a nitrogen content of not more than 0.002 % (mass fraction). Do not dry in an oven before use.

### Blank Tests

Carry out a blank test following the currently used procedure for digestion, distillation and titration taking 2 ml of water and about 0.7 g of sucrose instead of the test portion. Keep a record of blank values. If blank values change, identify the cause.

**Note: The amount of titrant used in the blank test should always be greater than 0.0 ml. Blanks within the same laboratory should be consistent across time.**

## 4.3 Recovery Tests

Regularly run recovery studies to check the accuracy of procedure and equipment:

- **Nitrogen loss** - Use 0.12 g ammonium sulphate and 0.67 g sucrose per flask weighted to the nearest 0.1 mg. Add all other reagents as stated in the method currently used (Kjeltabs,  $H_2SO_4$ , etc). Digest and distil under same conditions as for sample. Recoveries shall be >99 %.
- **Digestion efficiency** - Use a test portion of minimum 0.15 g of tryptophan or acetanilide and 0.67 g sucrose per flask weighted to the nearest 0.1 mg. Determine the nitrogen content according to the current procedure in use. The recoveries of tryptophan shall be >98.5 %; the recoveries of acetanilide shall be >99.5 %.
- **Distillation and titration efficiency** - Distil 0.10 - 0.15 g  $\pm 0.0001$  g ammonium sulphate, omitting the digestion step. The recoveries should be >99.5 %.

**Note: Results less than 98.5 % or more than 101.0 % in either of the recovery tests indicate failures in the procedure and/or inaccurate concentration of the standard volumetric hydrochloric acid solution (should be adjusted to four decimals accuracy according to procedure in AN 300)**

## Kjeltec™ 8100 Distillation Unit Tecator™ 2508/2520 Digestor

### 1 Scope

This PQ applies to the Digestion system 2508/2520 (including exhaust and scrubber unit) and Kjeltec 8100 Distillation Unit manufactured by FOSS Analytical. The user of the instrument performs the PQ.

### 2 Intended Use

The Digestion system (including exhaust and scrubber) and Kjeltec 8100 Distillation Unit are intended for laboratory use analyzing parameters as specified in FOSS Application Notes.

### 3 Purpose

The guidelines are intended to assist the user in successfully developing Performance Qualifications for the specific application(s) to which the instrument is applied.

The Performance Qualification (PQ) includes the process of demonstrating that the Digestion system 2508/2520 (including exhaust and scrubber unit) and the Kjeltec 8100 Distillation unit consistently perform according to a specification appropriate for its routine use. Main activities in the PQ phase are:

- Preventive maintenance
- On-going verification tests

This document suggests routines to fulfill the requirements for an acceptable PQ but the final procedure should be adapted to local routines for similar equipment.

### 4 Definition of Test Procedures

#### 4.1 Preventive Maintenance

Maintenance of the Kjeltec 8100 should be performed according to the instructions in manual, see User Manual Kjeltec 8100/8200 Distillation Unit, chapter 5. Maintenance. A yearly service is recommended (service agreement).

Maintenance of the Digestion block (including exhaust and scrubber) should be performed according to instruction in the user manual, see User Manual Tecator Digestor, chapter 5. Maintenance.

#### Dedicated Analytical Solutions

FOSS Analytical AS  
65 Slangerupgade  
DK-3400 Hillerød  
Denmark

Tel +45 7010 3370  
Fax +45 7010 3371  
E-mail support@foss.dk  
Web www.foss.dk

FOSS Analytical AB  
Box 70  
SE-263 21 Höganäs  
Sweden

Tel +46 42 361500  
Fax +46 42 360369  
E-mail support@foss.se  
Web www.foss.se

### External Quality Control Program

It is recommended to participate in an external quality control program, such a proficiency program or ring test, with equivalent sample material as analysed within the laboratory.

### Calculation and Expression of Results

$$W_N = \frac{14.007(V_s - V_b)N \times 100\%}{m}$$

Where:

$W_N$  is the nitrogen content of the sample, expressed as a percentage by mass.

$V_s$  is the numerical value of the volume of the hydrochloric acid standard volumetric solution used in the sample test, in milliliters, expressed to the nearest 0.05 ml.

$V_b$  is the numerical value of the volume of the hydrochloric acid standard volumetric solution used in the blank test, in milliliters, expressed to the nearest 0.05 ml.

$N$  is the numerical value of the exact normality of the hydrochloric acid standard volumetric solution, expressed to four decimal places.

$m$  is the numerical value of the mass of the test portion, in milligrams, expressed to the nearest 1 mg for sample weights >1 g or to the nearest 0.1 mg for sample weights <1 g.

## 5 Maintenance

### 5.1 Maintenance Kjeltec™ 8100

See instructions in User Manual - Kjeltec 8100/8200, chapter 5 Maintenance.

### 5.2 Maintenance Tecator™ Digestor

See instructions in User Manual - Tecator Digestor, chapter 5 Maintenance.

## 6 The Maintenance Record Charts

This record charts are provided to assist you in keeping your system in good working order. Please make copies and use them regularly as they can often help us to help you in the unlikely event a system malfunction.



### 6.2.5 Yearly Maintenance

[illegible]

## เอกสารไม่ควบคุม

### 6.2.3 Every 1-3 Months Maintenance

[illegible]

## เอกสารไม่ควบคุม

### 6.2.2.6 Exchange of Parts and Reagents Maintenance

[illegible]

## เอกสารไม่ควบคุม

#### 6.2.4 Additional Maintenance

[illegible]

## เอกสารไม่ควบคุม





Cert. No.: 23TM726  
Page : 1 of 3

## Certificate of Calibration

Equipment : Cooled Incubator  
Manufacturer : Binder  
Model : KB 400 E6  
Serial No. : 2020000015535  
ID No. : UAE.MIC.018/2564  
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260  
Location : Microbiology Laboratory (302)  
Received Order : 27 April 2023  
Calibration Date : 27 April 2023  
Ambient Temperature : ( 26 ± 10 ) °C  
Relative Humidity : ( 50 ± 30 ) %  
Calibrated by : Tawatchai Pama  
Approved by :   
( ) Pornthippa Tameyakul  
( ✓ ) Malee Butkrues  
( ) Suwit Imjai  
Issue Date : 12 May 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : Cooled Incubator  
Condition As-Received : Used Item  
Reference : 2304-0461OC-1  
Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector ( RTD ).  
The temperature scale used was based on ITS-90.

### Condition of this result of calibration

#### 1. Reference standard instrument-

Instrument	Model	Serial No.	Cert. No.	Due Date
1 ) Data Acquisition	34972A	MY57013711	22LM93	02 Jul 2023

2. This certificate is valid only to the item calibrated on date and place of calibration.


3. This certification is traceable to the International System of Unit.

Result of Calibration :- ( \* ) Without Adjustment

Function of UUC\* : Temperature Source

Fresh air setting : Close

Environment during calibration		
	Beginning	Finished
Temp. ( °C )	20	19
REL.Humid. ( % )	72	82
AC Supply ( Volt )	230	231

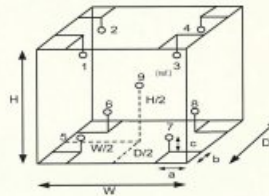

#### Probe Installation Details :

a = 10 cm  
b = 10 cm  
c = 10 cm

#### Dimension of Chamber :

D = 0.48 m  
W = 0.65 m  
H = 1.2 m  
Capacity = 0.37 m<sup>3</sup>

Position :	Ref. Std. ID No.:
1	22-18RTD-2/1
2	18RTD-2/2
3	18RTD-2/3
4	18RTD-2/4
5	18RTD-2/5
6	18RTD-2/6
7	18RTD-2/7
8	18RTD-2/8
9 (ref.)	18RTD-2/9

เอกสารไม่ควบคุม



Equipment : Cooled Incubator  
Condition As-Received : Used Item  
Reference : 2304-0461OC-1  
Result of Calibration :- ( \* ) Without Adjustment  
Function of UUC\* : Temperature Source  
Fresh air setting : Close

Cert. No.: 23TM726  
Page : 3 of 3

Calibration Point ( °C )	UUC* Setting ( °C )	UUC* Reading ( °C )	Temperature stability ( ± °C )	Temperature uniformity ( °C )	Overall Variation ( °C )	Coverage Factor k
35.0	35.0	35.0	0.0090	0.16	0.21	2

Calibration Point ( °C )	Measured Temperature ( °C )									Uncertainty ( ± °C )
	1	2	3	4	5	6	7	8	9 (ref.)	
35.0	34.913	34.997	34.834	34.893	35.034	35.027	35.025	35.035	34.980	0.30

Average\* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC\* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-000-

เอกสารไม่ควบคุม



Cert. No.: 23TM728  
Page : 1 of 3

## Certificate of Calibration

Equipment : Incubator  
Manufacturer : Memmert  
Model : IPP 260  
Serial No. : V616.0066  
ID No. : UAE.MIC.032/2559  
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260  
Location : Microbiology Laboratory (302)  
Received Order : 27 April 2023  
Calibration Date : 27 - 28 April 2023  
Ambient Temperature : ( 26 ± 10 ) °C  
Relative Humidity : ( 50 ± 30 ) %  
Calibrated by : Tawatchai Pama  
Approved by :   
( ) Pornthippa Tameyakul  
( ✓ ) Malee Butkrues  
( ) Suwit Imjai  
Issue Date : 11 May 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : Incubator  
Condition As-Received : Used Item  
Reference : 2304-0461OC-6  
Result of Calibration : ( \* ) Without Adjustment  
Function of UUC\* : Temperature Source  
Fresh air setting : Close

Cert. No.: 23TM728  
Page : 3 of 3

Calibration Point ( °C )	UUC* Setting ( °C )	UUC* Reading ( °C )	Temperature stability ( ± °C )	Temperature uniformity ( °C )	Overall Variation ( °C )	Coverage Factor k
25.0	25.0	25.0	0.020	0.81	1.2	2
36.0	36.0	36.0	0.15	1.1	1.6	2

Calibration Point ( °C )	Measured Temperature ( °C )									Uncertainty ( ± °C )
	1	2	3	4	5	6	7	8	9 (ref.)	
25.0	25.541	25.354	25.388	25.278	24.341	24.349	24.379	24.455	24.747	0.30
36.0	35.275	35.351	35.768	35.941	36.543	36.590	36.653	36.728	36.232	0.39

Average\* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC\* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-000-

เอกสารไม่ควบคุม

~ 4450000



Equipment : Incubator  
Condition As-Received : Used Item  
Reference : 2304-0461OC-6

Cert. No.: 23TM728  
Page : 2 of 3

Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector ( RTD ).

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1 ) Data Acquisition	34972A	MY57013711	22LM93	02 Jul 2023

2. This certificate is valid only to the item calibrated on date and place of calibration.


3. This certification is traceable to the International System of Unit.

Result of Calibration :- ( \* ) Without Adjustment

Function of UUC\* : Temperature Source

Fresh air setting : Close

Environment during calibration		
	Beginning	Finished
Temp. ( °C )	25	22
REL.Humid. ( % )	76	83
AC Supply ( Volt )	231	231

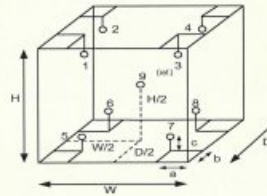

Position :	Ref. Std. ID No.:
1	22-18RTD-2/1
2	18RTD-2/2
3	18RTD-2/3
4	18RTD-2/4
5	18RTD-2/5
6	18RTD-2/6
7	18RTD-2/7
8	18RTD-2/8
9 (ref.)	18RTD-2/9

Probe Installation Details :

a = 10 cm  
b = 10 cm  
c = 10 cm

Dimension of Chamber :

D = 0.50 m  
W = 0.64 m  
H = 0.80 m  
Capacity = 0.26 m<sup>3</sup>

เอกสารไม่ควบคุม

~ 4450000



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES & EQUIPMENT CALIBRATION AND TESTING SERVICES  
534/4 PATTANAKARN ROAD SOI 16, SUANLIANG, SUANLIANG BANGKOK 10230  
TEL. 0-2717-3800-29 FAX. 0-2719-9484



Cert. No.: 23TM193  
Page : 1 of 3

## Certificate of Calibration

Equipment : Water Bath  
Manufacturer : Memmert  
Model : WNE 14  
Serial No. : L416.0606  
ID No. : UAE.MIC.002/2560  
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260  
Location : Microbiology Laboratory  
Received Order : 15 February 2023  
Calibration Date : 15 February 2023  
Ambient Temperature : ( 28 ± 10 ) °C  
Relative Humidity : ( 50 ± 30 ) %  
Calibrated by : Suwit Imjai  
Approved by :   
( ) Ponthippa Tameyakul  
( / ) Malee Butkruea

Issue Date : 24 February 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services & Equipment Calibration and Testing Services

เอกสารไม่ควบคุม



Equipment : Water Bath  
Condition As-Received : Used Item  
Reference : 2302-0295OC-2

Cert. No.: 23TM193  
Page : 2 of 3

Procedure Used :-

Calibration were conducted using in-house calibration procedure CP-OT04 according to direct measurement method with Data Acquisition which connected with Industrial Platinum Resistance Thermometer ( IPRT ).

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1 ) Data Acquisition	34972A	MY59003411	22LM165	26 Nov 2023


2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

Result of Calibration :- ( \* ) Without Adjustment

Function of UUC\* : Temperature Source

	Environmental		AC Voltage Supply
	( °C )	( %R.H. )	( Volt )
Beginning of Calibration	22	65	231
Finished of Calibration	23	61	231

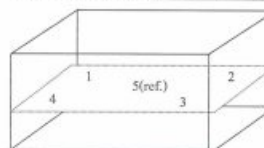

Front

Position :	Ref. Std. ID No.:
1	4804539-001
2	4804539-002
3	4804539-003
4	4804539-004
5 (ref.)	4804539-005

เอกสารไม่ควบคุม





Equipment : Water Bath  
Condition As-Received : Used Item  
Reference : 2302-0295OC-2  
Result of Calibration : ( \* ) Without Adjustment  
Function of UUC\* : Temperature Source

Cert. No.: 23TM193  
Page : 3 of 3

Calibration point ( °C )	UUC* Setting ( °C )	UUC* Reading ( °C )	Average* Standard Reading ( °C )				
			Position				
			1	2	3	4	5 (ref.)
44.5	44.5	44.5	44.453	44.437	44.428	44.477	44.459

Calibration point ( °C )	Uniformity ( °C )	Stability ( ± °C )	Uncertainty ( ± °C )	Coverage Factor k
44.5	0.079	0.036	0.15	2

Average\* : The average of 30 values in each position.

Uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Stability : One-half of the greatest maximum difference of measured temperature at any one probe.

UUC\* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity.

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor  $k$ , providing a level of confidence of approximately 95 %.

-000-

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES  
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250  
TEL. 0-2717-3000-29 FAX. 0-2719-9484



Cert. No.: 23TM194  
Page : 1 of 3

## Certificate of Calibration

Equipment : Water Bath  
Manufacturer : Memmert  
Model : WNE 14  
Serial No. : L416.0612  
ID No. : UAE.MIC.003/2560  
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260  
Location : Microbiology Laboratory  
Received Order : 15 February 2023  
Calibration Date : 15 February 2023  
Ambient Temperature : ( 26 ± 10 ) °C  
Relative Humidity : ( 50 ± 30 ) %  
Calibrated by : Suwit Imjai  
Approved by :   
( ) Pornthippa Tameyakul  
( ) Malee Butkruea  
Issue Date : 24 February 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services

เอกสารไม่ควบคุม



Equipment : Water Bath  
Condition As-Received : Used Item  
Reference : 2302-0295OC-3  
Procedure Used :-

Cert. No.: 23TM194  
Page : 2 of 3

Calibration were conducted using in-house calibration procedure CP-OT04 according to direct measurement method with Data Acquisition which connected with Industrial Platinum Resistance Thermometer (IPRT).

The temperature scale used was based on ITS-90.

### Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1 ) Data Acquisition	34972A	MY59003411	22LM165	26 Nov 2023


2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

Result of Calibration :- ( \* ) Without Adjustment

Function of UUC\* : Temperature Source

	Environmental		AC Voltage Supply
	( °C )	( %R.H. )	( Volt )
Beginning of Calibration	22	65	231
Finished of Calibration	22	63	230

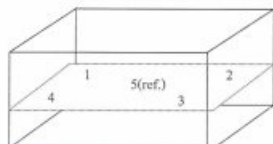

Front

Position :	Ref. Std. ID No.:
1	4804539-001
2	4804539-002
3	4804539-003
4	4804539-004
5(ref.)	4804539-005

เอกสารไม่ควบคุม



Equipment : Water Bath  
Condition As-Received : Used Item  
Reference : 2302-0295OC-3  
Result of Calibration : ( \* ) Without Adjustment  
Function of UUC\* : Temperature Source

Cert. No.: 23TM194  
Page : 3 of 3

Calibration point ( °C )	UUC* Setting ( °C )	UUC* Reading ( °C )	Average* Standard Reading ( °C )				
			Position				
			1	2	3	4	5 (ref.)
44.5	44.5	44.6	44.520	44.509	44.498	44.552	44.530

Calibration point ( °C )	Uniformity ( °C )	Stability ( ± °C )	Uncertainty ( ± °C )	Coverage Factor k
44.5	0.077	0.037	0.15	2

Average\* : The average of 30 values in each position.

Uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Stability : One-half of the greatest maximum difference of measured temperature at any one probe.

UUC\* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity.

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor  $k$ , providing a level of confidence of approximately 95 %.

-000-

เอกสารไม่ควบคุม





Cert. No.: 23TM763  
Page : 1 of 3

## Certificate of Calibration

Equipment : Autoclave  
Manufacturer : ALP  
Model : CL-40L  
Serial No. : 808763  
ID No. : UAE.MIC.026/2563  
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260  
Location : Microbiology Laboratory (301)  
Received Order : 27 April 2023  
Calibration Date : 27 April 2023  
Ambient Temperature : ( 26 ± 10 ) °C  
Relative Humidity : ( 50 ± 30 ) %  
Calibrated by : Preecha Hlahib  
Approved by :   
( ) Pornthippa Tameyakul  
( ) Melee Butkruea  
( ) Suwit Imjai  
Issue Date : 11 May 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3: Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม

A 0053944



Equipment : Autoclave  
Condition As-Received : Used Item  
Reference : 2304-0461OC-2  
Procedure Used :-

Cert. No.: 23TM763  
Page : 2 of 3

Calibration were conducted using in-house calibration procedure CP-OT03 according to direct measurement method with Data Acquisition which connected with Thermocouple Type T  
The temperature scale used was based on ITS-90.

### Condition of this result of calibration

#### 1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1 ) Data Acquisition	34972A	MY59003411	22LM165	26 Nov 2023

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.


4. This result of calibration covers laboratory autoclaves for the sterilization of goods and material which could be infected with organisms categorized as Hazard Group 1, 2 and 3\*\*

(\* = Categorization of pathogens according to hazard and categories of containment, second edition, 1990 )  
It does not cover autoclaves for use with material infect with organisms in Hazard Group 4, for which complete containment and sterilization of infected condensate is considered to be essential.

This result of calibration does not apply to sterilizers or disinfectors used for medical, dental, pharmaceutical or veterinary purposes which are directly concerned with patient care, or those used for fabrics subjected to sterilization which are required to be dry at the end of cycle.

Result of Calibration :- ( \* ) Without Adjustment

Function of UUC\* : Temperature Source

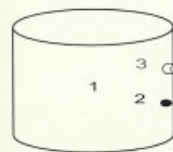

	Environmental		
	( °C )	( %R.H. )	( Volt )
Beginning of Calibration	27	60	220
Finished of Calibration	27	58	220

Position	Description	Ref. Std. ID No.:
1 =	Center of chamber	18-20TC-04
2 =	Temperature sensor	18-20TC-05
3 =	Exhaust port	18-20TC-06

เอกสารไม่ควบคุม

A 1159968



Equipment : Autoclave  
Condition As-Received : Used Item  
Reference : 2304-0461OC-2  
Result of Calibration :- ( \* ) Without Adjustment  
Function of UUC\* : Temperature Source

Cert. No.: 23TM763  
Page : 3 of 3

Operating parameter Set : Temperature = 115.0 °C  
Sterilization period = 15 minute

UUC* Setting ( °C )	UUC* Reading ( °C )	Position	Average* Standard Reading ( °C )	Stability ( ± °C )	Pressure Reading ( MPa )	Uncertainty ( ± °C )	Coverage Factor k
115.0	115.0	1	115.213	0.22	0.08	0.75	2
		2	115.168				
		3	115.260				

Operating parameter Set : Temperature = 121.0 °C  
Sterilization period = 30 minute

UUC* Setting ( °C )	UUC* Reading ( °C )	Position	Average* Standard Reading ( °C )	Stability ( ± °C )	Pressure Reading ( MPa )	Uncertainty ( ± °C )	Coverage Factor k
121.0	121.0	1	121.260	0.29	1.1	0.75	2
		2	121.224				
		3	121.284				

Average\* : The average of 30 values in each position.

Stability : One-half of the greatest maximum difference of measured temperature at any one probe.

UUC\* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-000-

เอกสารไม่ควบคุม

A 1159967



Cert. No.: 23TM763  
Page : 1 of 3

## Certificate of Calibration

Equipment : Autoclave  
Manufacturer : ALP  
Model : CL-40L  
Serial No. : 808763  
ID No. : UAE.MIC.026/2563  
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260  
Location : Microbiology Laboratory (301)  
Received Order : 27 April 2023  
Calibration Date : 27 April 2023  
Ambient Temperature : ( 26 ± 10 ) °C  
Relative Humidity : ( 50 ± 30 ) %  
Calibrated by : Preecha Hlahib  
Approved by :   
( ) Pornthippa Tameyakul  
( ) Melee Butkruea  
( ) Suwit Imjai  
Issue Date : 11 May 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3: Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม

A 0053944





Equipment : Autoclave  
Condition As-Received : Used Item  
Reference : 2304-0461OC-2  
Procedure Used :-

Cert. No.: 23TM763  
Page : 2 of 3

Calibration were conducted using in-house calibration procedure CP-OT03 according to direct measurement method with Data Acquisition which connected with Thermocouple Type T  
The temperature scale used was based on ITS-90.

#### Condition of this result of calibration

##### 1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1) Data Acquisition	34972A	MY59003411	22LM165	26 Nov 2023

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

4. This result of calibration covers laboratory autoclaves for the sterilization of goods and material which

could be infected with organisms categorized as Hazard Group 1, 2 and 3\*\*

(\*\* = Categorization of pathogens according to hazard and categories of containment, second edition, 1990)


It does not cover autoclaves for use with material infect with organisms in Hazard Group 4, for which

complete containment and sterilization of infected condensate is considered to be essential.

This result of calibration does not apply to sterilizers or disinfectors used for medical, dental, pharmaceutical or veterinary purposes which are directly concerned with patient care, or those used for fabrics subjected to sterilization which are required to be dry at the end of cycle.

Result of Calibration :- ( \* ) Without Adjustment

Function of UUC\* : Temperature Source

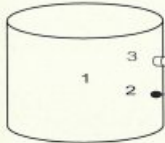

	Environmental		
	( °C )	( %R.H. )	( Volt )
Beginning of Calibration	27	60	220
Finished of Calibration	27	58	220

Position	Description	Ref. Std. ID No.
1 =	Center of chamber	18-20TC-04
2 =	Temperature sensor	18-20TC-05
3 =	Exhaust port	18-20TC-06

เอกสารไม่ควบคุม

a 1159968



Equipment : Autoclave  
Condition As-Received : Used Item  
Reference : 2304-0461OC-2

Cert. No.: 23TM763  
Page : 3 of 3

Result of Calibration :- ( \* ) Without Adjustment

Function of UUC\* : Temperature Source

Operating parameter Set : Temperature = 115.0 °C  
Sterilization period = 15 minute

UUC* Setting ( °C )	UUC* Reading ( °C )	Position	Average* Standard Reading ( °C )	Stability ( ± °C )	Pressure Reading ( MPa )	Uncertainty ( ± °C )	Coverage Factor k
115.0	115.0	1	115.213	0.22	0.08	0.75	2
		2	115.166				
		3	115.260				

Operating parameter Set : Temperature = 121.0 °C  
Sterilization period = 30 minute

UUC* Setting ( °C )	UUC* Reading ( °C )	Position	Average* Standard Reading ( °C )	Stability ( ± °C )	Pressure Reading ( MPa )	Uncertainty ( ± °C )	Coverage Factor k
121.0	121.0	1	121.260	0.29	1.1	0.75	2
		2	121.224				
		3	121.284				

Average\* : The average of 30 values in each position.

Stability : One-half of the greatest maximum difference of measured temperature at any one probe.

UUC\* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-000-

เอกสารไม่ควบคุม

a 1159967



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES & EQUIPMENT CALIBRATION AND TESTING SERVICES  
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250  
TEL. 0-2717-3000-29 FAX. 0-2719-9484



Cert. No.: 23TM376  
Page : 1 of 3

## Certificate of Calibration

Equipment : Refrigerator  
Manufacturer : Sanyo  
Model : SBC-337KD(GYN)  
Serial No. : 71100607  
ID No. : UAE.MIC.003/2551  
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10280  
Location : Microbiology Laboratory (302)  
Received Order : 11 April 2023  
Calibration Date : 12 April 2023  
Ambient Temperature : ( 26 ± 10 ) °C  
Relative Humidity : ( 50 ± 30 ) %  
Calibrated by : Preecha Hlahib  
Approved by :   
( ) Pornthippa Tameyakul  
( ) Maloe Butkrues  
( ) Suwit Imjai  
Issue Date : 24 April 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม

a 0057754



Equipment : Refrigerator  
Condition As-Received : Used Item  
Reference : 2304-0155OC-2  
Result of Calibration :- ( \* ) Without Adjustment  
Function of UUC\* : Temperature Source  
Fresh air setting : Not Available

Cert. No.: 23TM376  
Page : 3 of 3

Calibration Point ( °C )	UUC* Setting ( °C )	UUC* Reading ( °C )	Temperature stability ( ± °C )	Temperature uniformity ( °C )	Overall Variation ( °C )	Coverage Factor k
4	4	-	3.6	2.3	8.2	2

Calibration Point ( °C )	Measured Temperature ( °C )									Uncertainty ( ± °C )
	1	2	3	4	5	6	7	8	9 (ref.)	
4	5.544	4.470	4.383	4.041	3.949	3.838	3.457	3.471	3.891	4.3

Average\* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured

temperature at the reference location which are observed at the same time or at as close an observation time as

possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC\* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-000-

เอกสารไม่ควบคุม

a 1159770



Equipment : Refrigerator  
Condition As-Received : Used Item  
Reference : 2304-01550C-2  
Procedure Used :-

Cert. No.: 23TM376  
Page : 2 of 3

Calibration were conducted using calibration procedure CP-OT02 according to direct measurement method with Data Acquisition which connected with Thermocouple Type T.

The temperature scale used was based on ITS-90.

#### Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1) Data Acquisition	34972A	MY49001451	23LM27	25 Feb 2024

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.


**Result of Calibration :-** ( \* ) Without Adjustment

**Function of UUC\* :** Temperature Source

**Fresh air setting :** Not Available

Environment during calibration		
	Beginning	Finished
Temp. ( °C )	27	30
REL.Humid. ( % )	65	70
AC Supply ( Volt )	220	221

Position :	Ref. Std. ID No.:
1	20-19TC-01
2	20-19TC-02
3	20-19TC-03
4	20-19TC-04
5	20-19TC-05
6	20-19TC-06
7	20-19TC-07
8	20-19TC-08
9 (ref.)	20-19TC-09

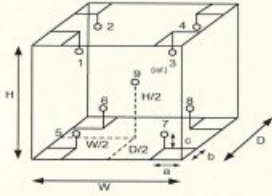

#### Probe Installation Details :

#### Dimension of Chamber :

a = 10 cm  
b = 10 cm  
c = 10 cm

D = 0.45 m  
W = 0.45 m  
H = 1.0 m  
Capacity = 0.20 m<sup>3</sup>

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES  
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250  
TEL. 0-2713-3080-27 FAX. 0-2719-9484



Cert. No.: 22TM1067  
Page.: 1 of 3

## Certificate of Calibration

Equipment : Refrigerator  
Manufacturer : Sanden Intercool  
Model : YPM11P  
Serial No. : YPM11P-150801805  
ID No. : UAE.MIC.012/2558  
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260  
Location : Microbiology Laboratory  
Received Order : 11 July 2022  
Calibration Date : 12 July 2022  
Ambient Temperature : ( 26 ± 10 ) °C  
Relative Humidity : ( 50 ± 30 ) %  
Calibrated by : Man Pattanapongsaiboon  
Approved by :   
( ) Pornthippa Tameysakul  
( ✓ ) Malee Butkruea  
( ) Suwit Imjai  
Issue Date : 18 July 2022

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written

Approval of the head of Corporate Services 3: Equipment Calibration and Testing

เอกสารไม่ควบคุม



Equipment : Refrigerator  
Condition As-Received : Used Item  
Reference : 2207-02450C-4  
Procedure Used :-

Cert. No.: 22TM1067  
Page.: 2 of 3

Calibration were conducted using calibration procedure CP-OT02 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector ( RTD ).

The temperature scale used was based on ITS-90.

#### Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1) Data Acquisition	34972A	MY57013823	22LM24	26 Feb 2023

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.


**Result of Calibration :-** ( \* ) Without Adjustment

**Function of UUC\* :** Temperature Source

**Fresh air setting :** Close

Environment during calibration		
	Beginning	Finished
Temp. ( °C )	25	25
REL.Humid. ( % )	56	59
AC Supply ( Volt )	222	221

Position :	Ref. Std. ID No.:
1	21-17RTD-01
2	21-17RTD-02
3	17RTD-03
4	17RTD-04
5	17RTD-05
6	17RTD-06
7	17RTD-07
8	17RTD-08
9 (ref.)	17RTD-09

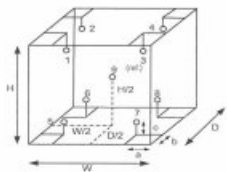

#### Probe Installation Details :

#### Dimension of Chamber :

a = 10 cm  
b = 10 cm  
c = 10 cm

D = 0.51 m  
W = 1.0 m  
H = 1.2 m  
Capacity = 0.61 m<sup>3</sup>

เอกสารไม่ควบคุม



Equipment : Refrigerator  
Condition As-Received : Used Item  
Reference : 2207-02450C-4  
Result of Calibration :- ( \* ) Without Adjustment  
Function of UUC\* : Temperature Source  
Fresh air setting : Close

Cert. No.: 22TM1067  
Page.: 3 of 3

Calibration Point ( °C )	UUC* Setting ( °C )	UUC* Reading ( °C )	Temperature stability ( ± °C )	Temperature uniformity ( °C )	Overall Variation ( °C )	Uncertainty ( ± °C )	Coverage Factor k
5.0	5.0	-1.4	2.5	3.3	5.6	3.1	2
Measured Temperature ( °C )							
Calibration Point ( °C )	Position						
5.0	1	2	3	4	5	6	7
	2.875	2.808	2.420	2.616	4.364	4.700	4.386
							4.129
							3.746

Average\* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor  
Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation

UUC\* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่ควบคุม





## Certificate of Calibration

Certificate No.: C01223732

Page: 2 of 2

Equipment: Balance  
 Model: PX623  
 Serial No. (or ID.): C236754745  
 Manufacturer: Ohaus  
 Condition: New

Certificate No.: C01223732  
 Issued Date: 09 December 2022  
 Job No.: KSPR2215576  
 Page: 1 of 2

Customer: United Analyst and Engineering Consultant Co., Ltd.  
 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak Sub-District,  
 Phrakhanong District, Bangkok, THAILAND 10260

Environment Condition: Temperature 26 °C ± 0.5 °C  
 Humidity 53 %RH ± 3.9 %RH

Calibration Place: United Analyst and Engineering Consultant Co., Ltd. (301 Microbiology Room)  
 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak Sub-District,  
 Phrakhanong District, Bangkok, THAILAND 10260

Calibration By: Mr. Adisai Maknoi  
 Calibration Date: 09 December 2022  
 The Method used: In-house method, CAL-WI-47, based on UKAS Lab 14  
 Traceability: This certificate is traceable to the SI Units maintained by National Institute of Metrology (NIMT), Thailand through DKSH Technology Co., Ltd. Certificate No. C02221765

## Calibration Results:

Without Adjustment

Eccentric Error: Weight to be 1/3 or 1/2 of Maximum capacity, taken from the center of the pan as a zero reference.

Nominal Test Value		200 (g)				
		Reference Points (g)				
		A	B	C	D	E
		-	0.000	0.000	0.000	0.000

Repeatability: Determination of the standard deviation of weighing balance., Readability 0.001 (g)

Nominal test value (g)	Standard Deviation
50	0.0004
500	0.0005

Error of Indication from nominal or conventional mass value., Readability 0.001 (g)

Nominal Value (g)	Conventional Mass (g)	Displayed Value (g)	Error of Indication (g)	Uncertainty (g)	k
1	1.0000	1.000	0.000	0.0010	2.03
5	5.0001	5.000	0.000	0.0010	2.03
10	10.0001	10.000	0.000	0.0010	2.03
20	20.0001	20.000	0.000	0.0010	2.03
50	50.0001	50.000	0.000	0.0010	2.03
100	100.0001	100.000	0.000	0.0011	2.03
200	200.0004	200.000	0.000	0.0011	2.02
300	300.0005	300.000	-0.001	0.0013	2.01
400	400.0008	400.001	0.000	0.0014	2.01
500	500.0003	500.000	0.000	0.0017	2.00
600	600.0004	600.000	0.000	0.0019	2.00

The End of Certificate

(Mr. Adisai Maknoi)  
 Person in charge

(Mr. Rungrod Jenkitrakulchai)  
 Authorized signatory

This certificate is issued the units of measurement according to the International System of Units (SI). It provides traceability of measurement to international or national standard or other recognized national standard laboratories.  
 The measurement uncertainty stated is the expanded uncertainty which is obtained from the standard uncertainty multiplied by the coverage factor (k=2) to provide a level of confidence of approximately 95%. It is determined in accordance with the Guide to Expression of Uncertainty in Measurement (GUM).  
 These results may be affected by deviations from specified conditions. The results relate only to the items tested, calibrated or sampled. The report shall not be reproduced except in full without approval of DKSH Technology Limited.

บริษัท ดีเคเอส อีเซีย จำกัด  
 DKSH Technology Limited  
 2533 ถนนสุขุมวิท แขวงคลองเตย เขตคลองเตย กรุงเทพมหานคร 10260  
 2533 Sukhumvit Road, Bangkok, Phrakhanong, Bangkok 10260  
 Phone: +66 2030 1100 Email: info.dksh@dksh.com Website: www.dksh.com/certificate-thailand

Delivering Growth - In Asia and Beyond.

CAL-FM-C01-14: 12 Sep 2022

เอกสารไม่ควบคุม

บริษัท ดีเคเอส อีเซีย จำกัด  
 DKSH Technology Limited  
 2533 ถนนสุขุมวิท แขวงคลองเตย เขตคลองเตย กรุงเทพมหานคร 10260  
 2533 Sukhumvit Road, Bangkok, Phrakhanong, Bangkok 10260  
 Phone: +66 2030 1100 Email: info.dksh@dksh.com Website: www.dksh.com/certificate-thailand

Delivering Growth - In Asia and Beyond.

CAL-FM-C01-14: 12 Sep 2022

เอกสารไม่ควบคุม

Refer to Certificate No.: C01223732

Page: 1 of 2

## Statements of conformity:

This conformity certificate documents the validity of the following statements of conformity based on the measurement results of corresponding calibration certificate:

The error of indication determined during calibration are under given measurement and environmental conditions and considering the expanded measurement uncertainty (coverage probability 95%) within the specification. The given measurement uncertainty already includes other all effects by according to the standard method, UKAS Lab14. Therefore, those parameters have not been assessed separately.

## Tolerance and Decision rules:

Assessment of the conformity of the measurement device are done based on direct comparison of the relevant measurement results with the tolerances and decision rule are prescribed by the customer.

- Decision rule:
- ☐ Choice A Binary Statement for Simple Acceptance Rule ( $w = 0$ ), Specific Risk < 50% PFA.
  - ☒ Choice B Non-binary statement with guard band ( $w = 1$  U), Pass or Fail Specific Risk < 2.5% PFA and Condition Pass or Condition Fail Specific Risk < 50% PFA.
  - ☐ Choice C Customer defined, Customers may define arbitrary multiple of  $r$  to have applied as guard band ( $w = r$  U).
- ; PFA - Probability of False Accept

(Mr. Rungrod Jenkitrakulchai)  
 Authorized signatory

Refer to Certificate No.: C01223732

Page: 2 of 2

## Statements of conformity:

Without Adjustment

Readability: 0.001 g

Nominal Value g	Error of Indication g	Guard band (w) g	Tolerance (±) g	Conformity
1	0.000	0.0010	0.002	Pass
5	0.000	0.0010	0.010	Pass
10	0.000	0.0010	0.020	Pass
20	0.000	0.0010	0.040	Pass
50	0.000	0.0010	0.100	Pass
100	0.000	0.0011	0.200	Pass
200	0.000	0.0011	0.400	Pass
300	-0.001	0.0013	0.600	Pass
400	0.000	0.0014	0.800	Pass
500	0.000	0.0017	1.000	Pass
600	0.000	0.0019	1.200	Pass

The validity of the statements of conformity cannot be guaranteed for different places of use, environmental conditions or improper use.

The End of Statements of conformity

บริษัท ดีเคเอส อีเซีย จำกัด  
 DKSH Technology Limited  
 2533 ถนนสุขุมวิท แขวงคลองเตย เขตคลองเตย กรุงเทพมหานคร 10260  
 2533 Sukhumvit Road, Bangkok, Phrakhanong, Bangkok 10260  
 Phone: +66 2030 1100 Email: info.dksh@dksh.com Website: www.dksh.com/certificate-thailand

Delivering Growth - In Asia and Beyond.

CAL-FM-C01-14: 12 Sep 2022


เอกสารไม่ควบคุม

บริษัท ดีเคเอส อีเซีย จำกัด  
 DKSH Technology Limited  
 2533 ถนนสุขุมวิท แขวงคลองเตย เขตคลองเตย กรุงเทพมหานคร 10260  
 2533 Sukhumvit Road, Bangkok, Phrakhanong, Bangkok 10260  
 Phone: +66 2030 1100 Email: info.dksh@dksh.com Website: www.dksh.com/certificate-thailand

Delivering Growth - In Asia and Beyond.

CAL-FM-C01-14: 12 Sep 2022

เอกสารไม่ควบคุม

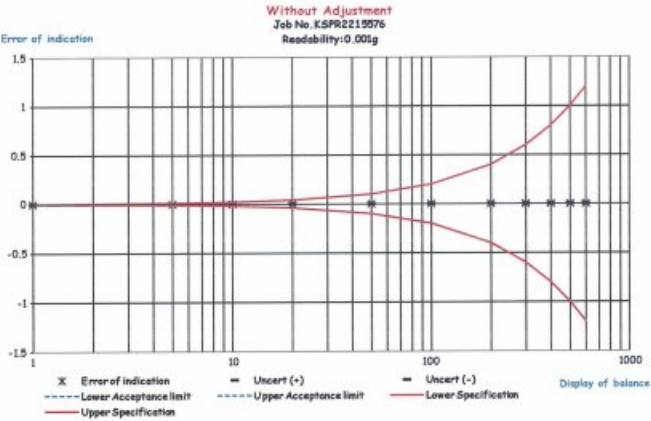

เอกสารไม่ควบคุม