

**ภาคผนวก จ**  
**เอกสารสอบเทียบเครื่องมือที่ใช้ในการตรวจวิเคราะห์**

รายการเครื่องมือหลักประจำห้องปฏิบัติการวิเคราะห์ สำหรับตรวจวิเคราะห์คุณภาพสิ่งแวดล้อม

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*	Remark
Laboratory Equipment for Air Analysis									
1	pH Meter	ความเป็นกรด-ด่าง	Mettler-Toledo	Seven Easy S20 / 1230525212	National Food Institute, Ministry of Industry, Thailand	2302181-001-01	24 Mar 23	22 Mar 24	-
2	pH Meter		Mettler-Toledo	SevenCompact S220/ C113432421	National Food Institute, Ministry of Industry, Thailand	2203527-001-01	5 Jul 22	4 Jul 23	-
3	Analytical Balance (Readability 0.01 mg)	ปริมาณของแข็งแขวนลอยทั้งหมด	Mettler-Toledo	XSR205DU / C210685394	Technology Promotion Association (Thailand-Japan)	23MM113	23 Sep 22	22 Sep 23	-
4	Hot Air Oven		Memmert	UF55 / B216.1666	Technology Promotion Association (Thailand-Japan)	22TM1490	19 Oct 22	18 Oct 23	-
5	Analytical Balance (Readability 0.1 mg)	น้ำมันและไขมัน	Mettler-Toledo	XSR204 / C117635043	National Food Institute, Ministry of Industry, Thailand	2202934-001-01	13 May 22	12 May 23	-
6	BOD Incubator		Arco	UC4-1320 / (UAE.WAO.015/2561)	Technology Promotion Association (Thailand-Japan)	23TM249	15 Feb 23	14 Feb 24	-
7	Incubator	ฟีดคอลโดลิทอรัมแบบเคีรีเย	Binder	BD 53 / 13-07343	Technology Promotion Association (Thailand-Japan)	23TM192	16 Feb 23	15 Feb 24	-
8	Water Bath		Memmert	WNE 14 / L416.0606	Technology Promotion Association (Thailand-Japan)	23TM193	15 Feb 23	14 Feb 24	-
9	Auto Clave		ALP	CL-40L / 807298	Technology Promotion Association (Thailand-Japan)	22TM1121	11 Jul 22	10 Jul 23	-
10	Analytical Balance		OHAUS	PX623 / C236754745	DKSH (Thailand) Ltd.	C01223732	9 Dec 22	8 Dec 23	-

# Calibration Certificate

Certificate No.: 2302181-001-01  
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.  
Address: 3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchack, Prakhnong, Bangkok 10260

Page 1 of 5

Equipment:	pH Meter
Manufacturer:	METTLER TOLEDO
Model:	SevenEasy pH
Serial No.:	1230525212
ID No.:	UAE.WAS.003/2553
Order No.:	2302181
Operation No.:	2302181-001
Date of Receipt:	14 March 2023
Date of Calibration:	24 March 2023

Calibrated by  
Mr. Phraphat Tuanjit  
Scientist  
Approved by  
Specialist, Division of Calibration Laboratory  
Responsible for the Technical Management Team  
Date of Issue: 24 March 2023

The uncertainties are for a confidence probability of approximately 95%.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65

©2008 บริษัทอสังหาริมทรัพย์ ไทยแลนด์ จำกัด (มหาชน) กรุงเทพมหานคร  
2008 Soi 35, Anuram Road, Bang Yi Khan Subdistrict, Bang Phlat District, Bangkok 10700, Thailand  
Tel : +66(0) 2422 8688 Fax : +66(0) 2422 8545



nfi orth

# Calibration Report

**Certificate No.:** 2302181-001-01  
**Equipment:** pH Meter  
**Manufacturer:**  
**Serial No.:**  
**ID No.:**

24 March 2023

**Location:** Chemical Calibration Laboratory, National Food Institute

Environment Condition: Ambient Temperature: (  $23.4 \pm 1.5$  ) °C Relative Humidity: (  $52 \pm 3$  ) %.

Condition of Equipment: Good Condition

### Condition of this Results of Calibration

1. Calibration Method  
In house method : W-CC-002 based on direct measurement by using standard voltage calibrator and certified reference material (CRM)

## 2. Reference Standards / Certified Reference Material

Instruments	Serial / ID No.	Manufacturer	Certificate No.	Due Date
2.1 DC Voltage Calibrator	2709007	Fluke	22E1959	17 June 2023
2.2 Digital Thermometer	2709007	Fluke	CC-650557-01	30 October 2023
2.3 Thermo-Hygro Meter	NILBH0347	PONPE	TE 650555-01	21 September 2023
Certified Reference Material	Lot No.	Manufacturer	Ref N	Expire Date
2.4 pH buffer 4.008 (Primary pH buffer Solution)	873608	CPAchem	PH216.L5	18 February 2023
2.5 pH buffer 6.885 (Primary pH buffer Solution)	873609	CPAchem	PH217.L5	18 February 2023
2.6 pH buffer 10.01 (Primary pH buffer Solution)	873611	CPAchem	PH220.L5	18 February 2023
2.7 pH buffer 7.00 (Standard pH buffer Solution)	873612	CPAchem	PH107.L5	18 February 2023

3. This certification is traceable to The International System of Unit (SI Unit)

3.1 Instruments No.2.1	through	NSC-TSI-TIS 17025 Laboratory Accreditation of Calibration No.0008
3.2 Instruments No.2.2	through	NSC-TSI-TIS 17025 Laboratory Accreditation of Calibration No.0061
3.3 Instruments No.2.3	through	NSC-TSI-TIS 17025 Laboratory Accreditation of Calibration No.0061
3.4 Certified Reference Material No. 2.4 to 2.6	traceable to	Primary measurement method- Harned cell using calibrated thermometer, barometer, and nanovoltmeter. The Standard Solution preparation and certified by CPAchem Ltd is accredited to ISO 17034 and ISO/IEC 17025
3.5 Certified Reference Material No 2.7	traceable to	BIM Ref'n HI-13 LoIn 25.05.2022; BIM Ref'n HI-16 LoIn 02.06.2022; BIM Ref'n HI-13 LoIn 25.05.2022; BIM Ref'n HI-16 LoIn 02.06.2022; the Standard Solution preparation and certified by CPAchem Ltd is accredited to ISO 17034 and ISO/IEC 17025

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

F-CS-012 Revision: 01 Date: 20-04-65

**เอกสารไมเคิล**

2008 บัณฑิตชั้นปริญญาโท สาขาวิชาศิลปกรรมศาสตร์  
2008 Sol 36, Anuramin Road, Bang Yi Khan Subdistrict, Bang Phlat District, Bangkok 10700, Thailand  
Tel : +66(0) 2422 8689 Fax : +66(0) 2422 8545



nfi orth

© 2008 บริษัท อีเอสเอส จำกัด 36 ถนนสุขุมวิท แขวงคลองตัน เขตคลองเตย กรุงเทพมหานคร 10110  
2008 Soi 36, Aun Anam Road, Bang Yi Khan Subdistrict, Bang Phai District, Bangkok 10700, Thailand  
Tel: +66(0) 2422 8688 Fax: +66(0) 2422 8545



of the





## Calibration Report

**Certificate No.:** 2302181-001-01  
**Equipment:** Digital Thermometer with RTD (pH Meter)  
Resolution: 0.1 °C Model: SevenEasy pH  
Serial No.: 1230525212 ID No.: UAE.WAS.003/2553  
Manufacturer: METTLER TOLEDO  
**Date of Calibration:** 24 March 2023

Page 5 of 5

**Calibration point:** 15.0, 25.0 and 30.0 °C  
**Calibration result:**

- The probe was immersed in liquid bath or dry bath to a minimum depth of 120 mm.
- Description of probe, model : N/A S/N : N/A
- Dimension of probe : Diameter 3 mm., Length 120 mm.,
- Sheath material : N/A

UUC* Reading (°C)	Standard Temperature (°C)	Correction Value (°C)	Uncertainty ± (°C)
15.2	14.999	- 0.2	0.12
25.2	24.999	- 0.2	0.12
30.2	29.999	- 0.2	0.12

### Note

- UUC\* : Unit Under Calibration

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k= 2, providing a level of confidence of approximately 95 %.

----- End -----

FCS-012 Revision: 01 Date: 20-04-65



เอกสารไม่ควบคุม  
nfi.or.th

## Calibration Certificate

**Certificate No.:** 2203527-001-01  
**Client name:** UNITED ANALYST AND ENGINEERING CONSULTANT CO., LTD.  
**Address:** 3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchack, Prakhong, Bangkok 10260

Page 1 of 5

**Equipment:** pH Meter  
**Manufacturer:** METTLER TOLEDO  
**Model:** Seven Compact S220  
**Serial No.:** C113432421

**ID No.:** UAE.WAT.009/2564

**Order No.:** 2203527

**Operation No.:** 2203527-001

**Date of Receipt:** 30 June 2022

**Date of Calibration:** 5 July 2022

**Calibrated by** Mr.Worapob Sooktong  
Scientist

**Approved by**

Manager, Division of Calibration Laboratory  
Responsible for the Technical Management Team

**Date of Issue:** 5 July 2022

The uncertainties are for a confidence probability of approximately 95%.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

FCS-009 Revision: 01 Date: 20-04-65



เอกสารไม่ควบคุม  
nfi.or.th





## Calibration Report

**Certificate No.:** 2203527-001-01  
**Equipment:** Digital Thermometer with RTD  
Resolution: 0.1 °C Model: Seven Compact S220  
Serial No.: C113432421 ID No.: UAE.WAT.009/2564  
Manufacturer: METTLER TOLEDO  
**Date of Calibration:** 5 July 2022

Page 4 of 5

**Location:** Calibration Laboratory, National Food Institute

**Environment Condition:** Ambient Temperature 25 °C ± 1 °C  
Relative Humidity 48 % ± 3 %

### Condition of this results of Calibration:

1. Calibration Method : - In house method: W-TE-025 by comparison with standard thermometer.  
- The Calibration is determined by comparing with a known temperature from a standard resistance thermometer.  
- The temperature scale in use at this laboratory is the International Temperature scale of 1990 ( ITS-90 ).
2. Reference Standard Instrument :

Instrument	Model	Serial No.	Certificate No.	Due Date	Through
HANDHELD THERMOMETER	1521	A85997	TE 650057-01	10-Dec-22	NATIONAL FOOD INSTITUTE
Platinum Resistance Thermometer (PRT)	385	509201			

Support Equipment : - Low Temperature Bath (ISOCAL-6), Model: Europa-6 Plus Basic, S/N: 341592/2

3. This certificate is traceable to International System of Units (SI Units).
4. This certificate was certified only for the instrument we calibrated.
5. This result of calibration was found accurate as shown on date and place of calibration only.
6. Condition of Calibrated item : Good
7. Result of Calibration : ☒ Without adjustment ☐ After adjustment

F-CS-012 Revision: 01 Date: 20-04-65

## Calibration Report

**Certificate No.:** 2203527-001-01  
**Equipment:** Digital Thermometer with RTD  
Resolution: 0.1 °C Model: Seven Compact S220  
Serial No.: C113432421 ID No.: UAE.WAT.009/2564  
Manufacturer: METTLER TOLEDO  
**Date of Calibration:** 5 July 2022

Page 5 of 5

**Calibration point:** 15.0, 25.0 and 35.0 °C

**Calibration result:** - The probe was immersed in liquid bath or dry bath to a minimum depth of 120 mm.  
- Description of probe, model : - S/N : -  
Dimension of probe : Diameter 9 mm., Length 120 mm.,  
Sheath material : Stainless Steel

UUC* Reading (°C)	Standard Temperature (°C)	Correction Value (°C)	Uncertainty ± (°C)
15.1	15.038	- 0.1	0.12
25.1	25.038	- 0.1	0.12
35.2	35.024	- 0.2	0.12

### Note

- UUC\* : Unit Under Calibration

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k= 2, providing a level of confidence of approximately 95 %.

----- End -----


F-CS-012 Revision: 01 Date: 20-04-65

Mettler-Toledo (Thailand) Ltd.  
846/4 - 846/5 Lasalle Rd., Bangna Tai Sub-District  
Bangna District, Bangkok 10260  
+66 2723 0382  
MT-TH.ServiceSupport@mt.com



## Accuracy Calibration Certificate

### Customer

Company:	United Analyst and Engineering Consultant Co., Ltd.		
Address:	3 Soi Udom Suk 41, Sukhumvit Rd., Bang Chak		
City:	Phra Khanong	Contact:	Suwit Chotnok
Zip / Postal:	10260		
State / Province:	Bangkok		
Order Number:	 * 0 3 3 2 5 3 8 1 3 *		

### Weighing Device

Manufacturer:	Mettler Toledo	Instrument Type:	Weighing Instrument
Model:	XPE205	Asset Number:	UAE.CAL.004/2561
Serial No.:	B748058497	Terminal Model:	PEAT
Building:	N/A	Terminal Serial No.:	B748058497
Floor:	4	Terminal Asset No.:	N/A
Room:	Laboratory Calibrate Glassware		

Range	Max. Capacity	Readability (d)
1	220 g	0.00001 g

### Procedure



Calibration Guideline:  
METTLER TOLEDO Work Instruction:  
EURAMET cg-18 v. 4.0 (11/2015)  
CP/W002/20

This calibration certificate contains measurements for As Found calibration. No As Left calibration was performed because the device was not modified after As Found calibration. Therefore, results for As Left correspond to As Found.

The sensitivity/span of the weighing instrument was adjusted before calibration with a built-in weight.


In accordance with EURAMET cg-18 (11/2015), the test loads were selected to reflect the specific use of the weighing device or to accommodate specific calibration conditions.

As Found	Temperature		Humidity	
	Start: 22.8 °C	End: 22.5 °C	Start: 62.2 %	End: 58.6 %


As Found Calibration Date:	23-Sep-2022	Calibrator:	
As Left Calibration Date:	N/A		
Issue Date:	26-Sep-2022	Approved Signatory:	 Suwicha Choykamchu
			Technical Manager / Head of Calibration Center

## Measurement Results

### Repeatability

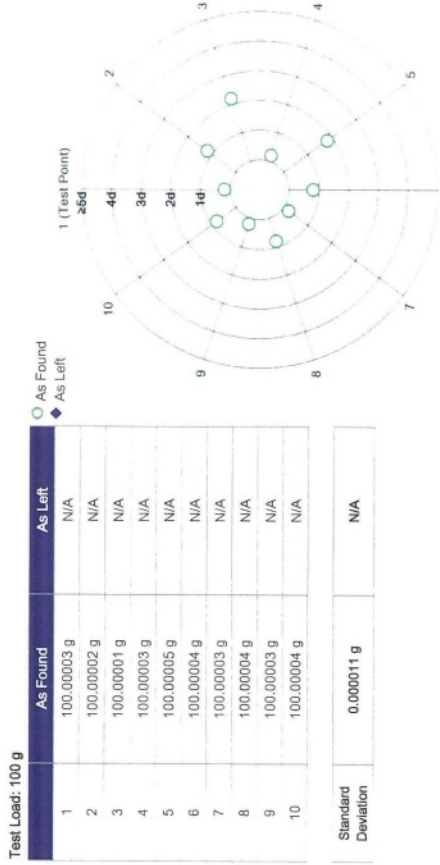

### Eccentricity

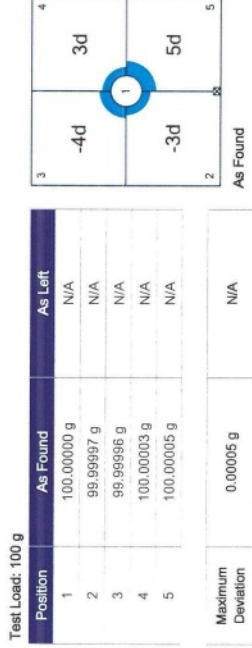




Error of Indication

As Found	Reference Value	Indication	Error of Indication	Expanded Uncertainty	k
1	0.00000 g	0.00000 g	0.00000 g	0.023 mg	2
2	9.99997 g	9.99997 g	0.00000 g	0.063 mg	2
3	29.99998 g	29.99999 g	0.00001 g	0.14 mg	2
4 <sup>*</sup>	49.99998 g	50.00000 g	0.00002 g	0.12 mg	2
5 <sup>*</sup>	69.99999 g	70.00000 g	0.00001 g	0.20 mg	2
6	90.00002 g	90.00002 g	0.00000 g	0.27 mg	2
7 <sup>*</sup>	100.00002 g	100.00005 g	0.00003 g	0.20 mg	2
8 <sup>*</sup>	120.00003 g	120.00005 g	0.00002 g	0.31 mg	2
9 <sup>*</sup>	150.00001 g	150.00002 g	0.00001 g	0.31 mg	2
10	170.00001 g	170.00007 g	0.00006 g	0.38 mg	2
11	200.00002 g	200.00000 g	-0.00002 g	0.35 mg	2

<sup>\*</sup>The calculated uncertainty was replaced by the CMC (Calibration and Measurement Capabilities) value because the calculated uncertainty was smaller than the CMC value.

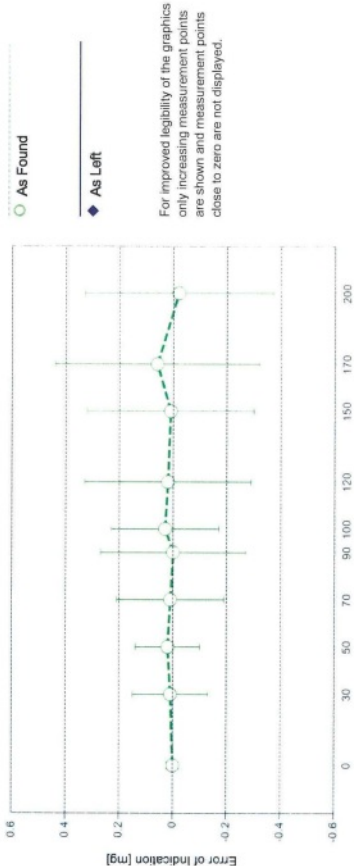

For improved legibility of the graphics only increasing measurement points are shown and measurement points close to zero are not displayed.

The uncertainty stated is the expanded uncertainty at calibration obtained by multiplying the standard combined uncertainty by the coverage factor k – which can be larger than 2 according to EURAMET cg-18. The value of the measurand lies within the assigned range of values with a probability of approximately 95%.

The user is responsible for maintaining environmental conditions and the settings of the weighing instrument when it was calibrated.

Test Equipment

All weights used for metrological testing are traceable to national or international standards. The weights were calibrated and certified by an accredited calibration laboratory.

Weight Set 1: OIML E2

Weight Set No.:	WS03	Date of Issue:	21-Sep-2021
Certificate Number:	175498	Calibration Due Date:	14-Mar-2023

Weight Set 2: OIML E2

Weight Set No.:	WS66	Date of Issue:	21-Oct-2021
Certificate Number:	C142784709	Calibration Due Date:	17-Apr-2023

Weight Set 3: OIML E2

Weight Set No.:	WS70	Date of Issue:	21-Oct-2021
Certificate Number:	C142784702	Calibration Due Date:	19-Mar-2023

Thermo Hygrometer

Equipment No.:	IN281	Date of Issue:	23-May-2022
Certificate Number:	22H1057	Calibration Due Date:	15-May-2023

Remarks

- FACT adjustment functionality activated
- Equipment condition: Good
- Next calibration according to customer's procedure
- Calibration data not decide by calibration laboratory.

End of Accredited Section

The information below and any attachments to this calibration certificate are not part of the accredited calibration.

Measurement Uncertainty of the Weighing Instrument in Use

Stated is the expanded uncertainty with  $k=2$  in use. The formula shall be used for the estimation of the uncertainty under consideration of the errors of indication. The value R represents the net load indication in the unit of measure of the device.

Temperature coefficient for the evaluation of the measurement uncertainty in use:  $1.0 \cdot 10^{-6} / K$

Temperature range on site for the evaluation of the measurement uncertainty in use: 3 K



Linearization of Uncertainty Equation

1	Range		As Found	As Left
	d	Max		
1	0.00001 g	220 g	$U_1 = 0.024 \text{ mg} + 0.00509 \text{ mg/g} \cdot R$	N/A

To optimize the stability of the linearization, besides of the zero load only increasing measurement points with a test load of 5% of the measurement range or larger are taken for the calculation of the linear equation.

Absolute and Relative Measurement Uncertainty in Use for Various Net Indications (Examples)

Net Indication	As Found	As Left
0.00220 g	0.024 mg	1.1%
0.02200 g	0.024 mg	0.11%
0.22000 g	0.025 mg	0.011%
2.20000 g	0.035 mg	0.0016%
220.00000 g	1.1 mg	0.00052%

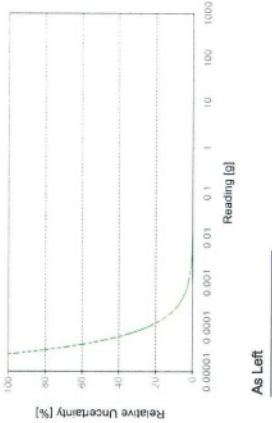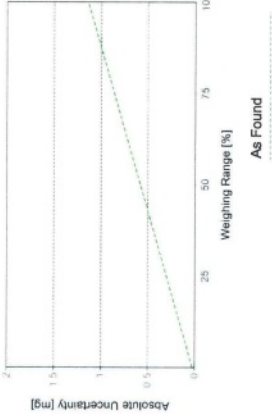

As Found

As Left



AS Found



AS Left



The weighing device meets the given process requirements.

The weighing device meets the given process requirements.

Tests Performed:



As Found



As Left

No adjustments/modifications made, As Left results correspond to As Found.


Process Requirements

Weighing Tolerance: 0.5%

Smallest Net Weight: 10.00000 g

Safety Factor: 2

Safe Weighing Range

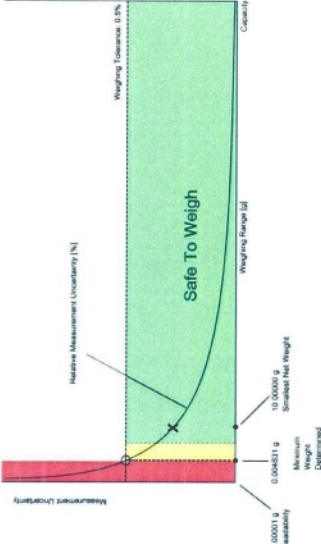

While the values in this graph reflect the actual calibration results, the measurement uncertainty curves are simply a visual representation. This graph reflects As Left testing, unless only As Found was performed.

Minimum Weight

As Found Minimum Weight Table

Tolerance	Minimum weights for different weighing tolerances and safety factors				
	Safety Factor				
	1	2	3	5	10
0.1%	0.024253 g	0.048754 g	0.073509 g	0.123795 g	0.254223 g
0.2%	0.012095 g	0.024253 g	0.036472 g	0.061100 g	0.123795 g
0.5%	0.004831 g	0.009671 g	0.014522 g	0.024253 g	0.048754 g
1%	0.002414 g	0.004831 g	0.007250 g	0.012095 g	0.024253 g
2%	0.001207 g	0.002414 g	0.003622 g	0.006040 g	0.012095 g
5%	0.000483 g	0.000965 g	0.001448 g	0.002414 g	0.004831 g

Pass: The determined minimum weight meets the requirement for the smallest net weight.

As Left Minimum Weight Table

Tolerance	Minimum weights for different weighing tolerances and safety factors				
	Safety Factor				
	1	2	3	5	10
0.1%	0.024253 g	0.048754 g	0.073509 g	0.123795 g	0.254223 g
0.2%	0.012095 g	0.024253 g	0.036472 g	0.061100 g	0.123795 g
0.5%	0.004831 g	0.009671 g	0.014522 g	0.024253 g	0.048754 g
1%	0.002414 g	0.004831 g	0.007250 g	0.012095 g	0.024253 g
2%	0.001207 g	0.002414 g	0.003622 g	0.006040 g	0.012095 g
5%	0.000483 g	0.000965 g	0.001448 g	0.002414 g	0.004831 g

Pass: The determined minimum weight meets the requirement for the smallest net weight.

At these net minimum weight values, the measurement uncertainty of the weighing device is equal to or less than 1/1 (no safety factor), 1/2, 1/3, 1/5, or 1/10 of the required tolerance. The values are calculated with  $k = 2$  and based on the linear formula of the measurement uncertainty of the weighing device in use.

The safety factor for As Found is always 1. This implies no safety factor. As Found testing looks at the behavior of the instrument from the past until test occurred. For the past, it is necessary to know that the tolerance was met, but not the safety factor. The safety factor is a proactive measure to apply for future measurements.

Notes on minimum weight values in above table:

1. If "N/A" is shown above, no appropriate value could be calculated.
2. METTLER TOLEDO is not responsible for the definition of the process requirements.

Measurement Results

Results Summary

	Repeatability		Eccentricity	Error of Indication
	As Found	As Left		
	✓	✓	✓	✓
	✓	✓	✓	✓

✓ = Passed  
✗ = Failed  
⚠ = Safety Factor not met

Repeatability

Test Load: 100 g

Tolerance	Control Limit	As Found		As Left	
		Std. Deviation	Result	Std. Deviation	Result
0.1%	0.005000 g		✓		✓
0.2%	0.010000 g		✓		✓
0.5%	0.025000 g	0.000011 g	✓	0.000011 g	✓
1%	0.050000 g		✓		✓
2%	0.100000 g		✓		✓
5%	0.250000 g		✓		✓

The weighing tolerance is met if the standard deviation is less than or equal to the corresponding control limit.

Eccentricity

Test Load: 100 g

Tolerance	Control Limit	As Found		As Left	
		Deviation	Result	Deviation	Result
0.1%	0.05000 g		✓		✓
0.2%	0.10000 g		✓		✓
0.5%	0.25000 g	0.00005 g	✓	0.00005 g	✓
1%	0.50000 g		✓		✓
2%	1.00000 g		✓		✓
5%	2.50000 g		✓		✓

The weighing tolerance is met if the deviation is less than or equal to the corresponding control limit.



Error of Indication

As Found

Control limits for various weighing tolerances									
Reference Value	Error	0.1%	0.2%	0.5%	1%	2%	5%		
0.00000 g	0.00000 g	N/A	N/A	N/A	N/A	N/A	N/A		
29.99998 g	0.00001 g	0.01500 g	0.03000 g	0.07500 g	0.15000 g	0.30000 g	0.75000 g		
49.99998 g	0.00002 g	0.02500 g	0.05000 g	0.12500 g	0.25000 g	0.50000 g	1.25000 g		
69.99999 g	0.00001 g	0.03500 g	0.07000 g	0.17500 g	0.35000 g	0.70000 g	1.75000 g		
90.00002 g	0.00000 g	0.04500 g	0.09000 g	0.22500 g	0.45000 g	0.90000 g	2.25000 g		
100.00002 g	0.00003 g	0.05000 g	0.10000 g	0.25000 g	0.50000 g	1.00000 g	2.50000 g		
120.00003 g	0.00002 g	0.06000 g	0.12000 g	0.30000 g	0.60000 g	1.20000 g	3.00000 g		
150.00001 g	0.00001 g	0.07500 g	0.15000 g	0.37500 g	0.75000 g	1.50000 g	3.75000 g		
170.00001 g	0.00006 g	0.08500 g	0.17000 g	0.42500 g	0.85000 g	1.70000 g	4.25000 g		
200.00002 g	-0.00002 g	0.10000 g	0.20000 g	0.50000 g	1.00000 g	2.00000 g	5.00000 g		
Result		✓	✓	✓	✓	✓	✓		✓

As Left

Control limits for various weighing tolerances									
Reference Value	Error	0.1%	0.2%	0.5%	1%	2%	5%		
0.00000 g	0.00000 g	N/A	N/A	N/A	N/A	N/A	N/A		
29.99998 g	0.00001 g	0.01500 g	0.03000 g	0.07500 g	0.15000 g	0.30000 g	0.75000 g		
49.99998 g	0.00002 g	0.02500 g	0.05000 g	0.12500 g	0.25000 g	0.50000 g	1.25000 g		
69.99999 g	0.00001 g	0.03500 g	0.07000 g	0.17500 g	0.35000 g	0.70000 g	1.75000 g		
90.00002 g	0.00000 g	0.04500 g	0.09000 g	0.22500 g	0.45000 g	0.90000 g	2.25000 g		
100.00002 g	0.00003 g	0.05000 g	0.10000 g	0.25000 g	0.50000 g	1.00000 g	2.50000 g		
120.00003 g	0.00002 g	0.06000 g	0.12000 g	0.30000 g	0.60000 g	1.20000 g	3.00000 g		
150.00001 g	0.00001 g	0.07500 g	0.15000 g	0.37500 g	0.75000 g	1.50000 g	3.75000 g		
170.00001 g	0.00006 g	0.08500 g	0.17000 g	0.42500 g	0.85000 g	1.70000 g	4.25000 g		
200.00002 g	-0.00002 g	0.10000 g	0.20000 g	0.50000 g	1.00000 g	2.00000 g	5.00000 g		
Result		✓	✓	✓	✓	✓	✓		✓

The weighing tolerance is met if the error (of indication) for each test point is less than or equal to the corresponding control limit for that particular weighing tolerance. Results at or close to the zero point cannot be assessed.



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES  
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250  
TEL. 0-2717-3000-27 FAX. 0-2719-4484



NSC-TB-17817025  
CALIBRATION 0008

Cert. No.: 22TM1490  
Page : 1 of 3

## Certificate of Calibration

Equipment : Hot Air Oven

Manufacturer : Memmert

Model : UF 55

Serial No. : B216.1666

ID No. : UAE.WAO.027/2559

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260  
Lab Floor 2

Received Order : 19 October 2022  
Calibration Date : 19 October 2022  
Ambient Temperature : ( 26 ± 10 ) °C  
Relative Humidity : ( 50 ± 30 ) %

Calibrated by : Preecha Hlahib

Approved by :   
Approved Signatory

( ) Pornthippa Tameyakul  
( ) Malee Butkruea  
(✓) Suwit Imjai

Issue Date : 31 October 2022

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.





Equipment : Hot Air Oven  
Condition As-Received : Used Item  
Reference : 2210-0575OC-1

Cert. No.: 22TM1490  
Page : 2 of 3

#### Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD ) and Thermocouple Type T.

The temperature scale used was based on ITS-90.

#### Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1) Data Acquisition	34970A	MY41021843	22LM4	10 Jan 2023


2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

Result of Calibration :- ( \* ) Without Adjustment

Function of UUC\* : Temperature Source

Fresh air setting : Close

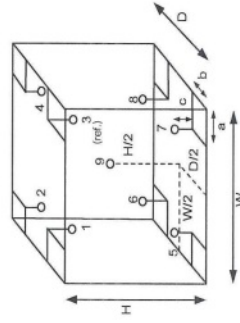

#### Probe Installation Details :

a =	5.0	cm
b =	5.0	cm
c =	5.0	cm
D =	0.33	m
W =	0.40	m
H =	0.40	m
Capacity =	0.053	m <sup>3</sup>

#### Dimension of Chamber :

Environment during calibration		
	Beginning	Finished
Temp. ( °C )	29	30
REL.Humid. ( % )	47	40
AC Supply ( Volt )	221	220

Ref. Std. ID No.: @ Calibration Point		
Position :	( 104 ) °C	( 140,180 ) °C
1	18-04RTD-01	21-04TC-01
2	18-04RTD-02	21-04TC-02
3	18-04RTD-03	21-04TC-03
4	18-04RTD-04	21-04TC-04
5	18-04RTD-05	21-04TC-05
6	18-04RTD-06	21-04TC-06
7	18-04RTD-07	21-04TC-07
8	18-04RTD-08	21-04TC-08
9 (ref.)	18-04RTD-09	21-04TC-09



Equipment : Hot Air Oven  
Condition As-Received : Used Item  
Reference : 2210-0575OC-1  
Result of Calibration :- ( \* ) Without Adjustment  
Function of UUC\* : Temperature Source  
Fresh air setting : Close

Cert. No.: 22TM1490  
Page : 3 of 3

Calibration Point ( °C )	UUC* Setting ( °C )	UUC* Reading ( °C )	Temperature stability ( ± °C )	Temperature uniformity ( °C )	Overall Variation ( °C )	Uncertainty ( ± °C )	Coverage Factor k
104.0	104.0	104.0	0.061	1.3	1.7	0.42	2
140.0	140.0	140.0	0.14	2.3	2.4	1.1	2
180.0	180.0	180.0	0.21	3.5	3.6	1.3	2

Measured Temperature ( °C )								
Position								
1	2	3	4	5	6	7	8	9 ( ref. )
103.076	103.876	103.777	104.124	104.667	104.426	104.012	103.928	104.370
138.199	139.189	138.808	139.550	140.266	139.622	139.293	139.385	140.369
177.930	179.267	178.643	179.753	181.011	180.093	179.496	179.743	181.278

Average\* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor  
Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation  
UUC\* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-o0o-



Equipment : Hot Air Oven  
Condition As-Received : Used Item  
Reference : 2210-0575OC-1

Cert. No.: 22TM1490  
Page : 2 of 3

#### Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD ) and Thermocouple Type T.

The temperature scale used was based on ITS-90.

#### Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1) Data Acquisition	34970A	MY41021843	22LM4	10 Jan 2023


2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

Result of Calibration :- ( \* ) Without Adjustment

Function of UUC\* : Temperature Source

Fresh air setting : Close

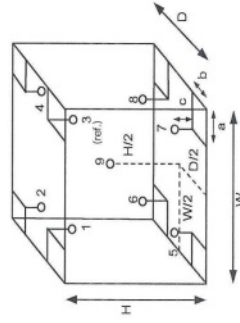

#### Probe Installation Details :

a =	5.0	cm
b =	5.0	cm
c =	5.0	cm
D =	0.33	m
W =	0.40	m
H =	0.40	m
Capacity =	0.053	m <sup>3</sup>

#### Dimension of Chamber :

Environment during calibration		
	Beginning	Finished
Temp. ( °C )	29	30
REL.Humid. ( % )	47	40
AC Supply ( Volt )	221	220

Ref. Std. ID No.: @ Calibration Point		
Position :	( 104 ) °C	( 140,180 ) °C
1	18-04RTD-01	21-04TC-01
2	18-04RTD-02	21-04TC-02
3	18-04RTD-03	21-04TC-03
4	18-04RTD-04	21-04TC-04
5	18-04RTD-05	21-04TC-05
6	18-04RTD-06	21-04TC-06
7	18-04RTD-07	21-04TC-07
8	18-04RTD-08	21-04TC-08
9 (ref.)	18-04RTD-09	21-04TC-09

เอกสารไม่ควบคุม

a 1133252

เอกสารไม่ควบคุม

a 1133251

## Calibration Certificate

**Certificate No.:** 2202934-001-01  
**Client name:** UNITED ANALYST AND ENGINEERING CONSULTANT CO., LTD.  
**Address:** 3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchack, Prakhnong, Bangkok 10260

Page 1 of 4

**Equipment:** Electronic Balance  
**Manufacturer:** METTLER TOLEDO  
**Model:** XSR204  
**Serial No.:** C117635043  
**ID No.:** UAE.WAS.012/2564  
**Order No.:** 2202934  
**Operation No.:** 2202934-001  
**Date of Receipt:** 13 May 2022  
**Date of Calibration:** 13 May 2022

**Calibrated by** Mr. Manas Somsak  
Specialist  
**Date of Issue:** 25 May 2022  
**Approved by** (Mr. [REDACTED])  
Manager, Division of Calibration Laboratory  
Responsible for the Technical Management Team

**The uncertainties are for a confidence probability of approximately 95%**

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

FCS-009 Revision: 01 Date: 20-04-65

## Calibration Report

**Certificate No.:** 2202934-001-01  
**Equipment:** Electronic Balance  
**Manufacturer:** METTLER TOLEDO  
**Model:** XSR204  
**Resolution:** 0.0001 g  
**Serial No.:** C117635043  
**ID No.:** UAE.WAS.012/2564  
**Capacity:** 220 g

Page 2 of 4

**Date of Calibration:** 13 May 2022  
**Environment Condition:** Ambient Temperature 22.3 ± 0.1 °C Relative Humidity: 47 ± 3 %  
**Place of Calibration:** Balance room (Water Analysis Unit), UNITED ANALYST AND ENGINEERING CONSULTANT CO., LTD.  
**Condition of Equipment:** Good Condition

### Condition of This Results of Calibration:

1. Calibration Method: NFI Method W-MA-001 In-House Method based on UKAS Lab 14 : 2019

2. Reference Standards:

Reference Standard	Model	Serial No	Calibrated By	Certificate No.	Due Date
Standard Weight Class E2	1mg to 200g	B505567572	TCS	M22041375	23 April 2023
Instrument	Model	Serial No	Calibrated By	Certificate No.	Due Date
Thermo-Hygro Meter	PONPE 490	NFLBTH 010/18	Quality Reborn	QR22-0350	18 February 2023

3. This certification is traceable to SI UNIT

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

### Calibration Results:




#### 1. Repeatability of Reading:

Nominal Value (g)	Standard Deviation of Reading (g)
100	0.000033
200	0.000032

#### 2. Off-Center Error:

A mass of 50 g was placed and moved to various position on pan.

The balance reading obtained is given in the table.

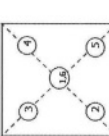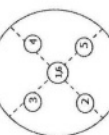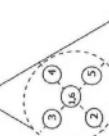						
<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>				
1 ( g )	2 ( g )	3 ( g )	4 ( g )	5 ( g )	6 ( g )	(Maximum Difference) ( g )
50 000	50 000	50 000	50 000	50 000	50 000	0 000

FCS-012 Revision: 01 Date: 20-04-65



# Calibration Report

Certificate No.: 2202934-001-01

**Equipment:**

### Electronic Balance

**Manufacturer:** METTLER TOLEDO

Model: XSR204

Resolution: 0.0001 g

Serial No.: C117635043

ID No.: UAE-WAS.012/2564

Capacity: 220 g

Date of Calibration: 13 May 2022

Page 3 of 4

**Calibration Results:** (Continued)

Calibration Results: (Continued)

**Calibration Range:** 0 - 200 g

**Calibration Adjustment: Internal Calibration**

3. Departure from Nominal Value:

Nominal Value ( g )	Standard Value ( g )	Average Reading ( g )	Correction ( g )	Uncertainty ( ± g )	Coverage Factor k
Unload	0.00000	0.0000	0.0000	0.000085	2.00
0.01	0.01000	0.0100	0.0000	0.000085	2.00
0.02	0.02000	0.0200	0.0000	0.000085	2.00
0.05	0.05000	0.0500	0.0000	0.000085	2.00
0.1	0.10001	0.1000	0.0000	0.000085	2.00
0.2	0.20001	0.2000	0.0000	0.000085	2.00
0.5	0.50002	0.5000	0.0000	0.000085	2.00
1	1.00001	1.0000	0.0000	0.000086	2.00
2	2.00003	2.0000	0.0000	0.000086	2.00
3	3.00004	3.0000	0.0000	0.000087	2.00
5	5.00002	5.0000	0.0000	0.000087	2.00
10	10.00001	10.0000	0.0000	0.000088	2.00
20	20.00004	20.0000	0.0000	0.000092	2.00
30	30.00005	30.0001	-0.0001	0.00010	2.00
40	40.00008	40.0001	0.0000	0.00011	2.00
45	45.00010	45.0001	0.0000	0.00013	2.00

F-CS-012 Revision: 01 Date: 20-04-65

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor  $k$ , providing a level of confidence of approximately 95 %.

----- End -----

F-CS-012 Revision: 01 Date: 20-04-65

2008 บะเข็สนวนสินธุ์ 35 นามะพรหมสินธุ์ แขวงท่าคันโท ขอนแก่น ๔๒1 กรุงเทพมหานคร  
 2008 Sol 35 Avun Ananin Road, Bang Yi Khan Subdistrict, Bang Phat District, Bangkok 10700, Thailand  
 Tel: +66(0) 2422 8568 Fax: +66(0) 2422 8545

2008 ขอเชิญคุณมาร่วมงานนี้! **เอกสารไม่ควบคุม**

2008 Soi 36, Arun Amarin Road, Bang Yi Khan Subdistrict, Bang Phlat District, Bangkok 10700, Thailand

Tel: +66(0) 2422 8688 Fax: +66(0) 2422 8545

nlorth



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES  
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250  
TEL: 0-2717-3000-29 FAX: 0-2719-9484



NSC-TSI-TS17025  
CALIBRATION 0008

Cert. No.: 23TM249  
Page : 1 of 3

## Certificate of Calibration

Equipment : BOD Incubator

Manufacturer : Arco

Model : UC4-1320

Serial No. : 13URC4S013201

ID No. : UAE.WAO.015/2561

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260  
Lab Floor 2

Received Order : 15 February 2023  
Calibration Date : 15 February 2023  
Ambient Temperature :  $(26 \pm 10) ^\circ\text{C}$   
Relative Humidity :  $(50 \pm 30) \%$

Calibrated by : Preecha Hlahib

Approved by :   
Ponthippa Tameyakul  
Malee Bulkruea  
Suwit Imjai

Issue Date : 24 February 2023

The Uncertainties are for a confidence probability of approximately 95 %

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม

A 0051476



Equipment : BOD Incubator  
Condition As-Received : Used Item  
Reference : 2302-0297OC-1  
Cert. No.: 23TM249  
Page : 2 of 3

### Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD ).

The temperature scale used was based on ITS-90.

### Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1 ) Data Acquisition	34972A	MYS7013711	22LM93	02 Jul 2023

2. This certificate is valid only to the item calibrated on date and place of calibration.


3. This certification is traceable to the International System of Unit.

Result of Calibration :- ( \* ) Without Adjustment

Function of UUC\* : Temperature Source

Fresh air setting : Not Available

Environment during calibration	
Beginning	Finished
Temp. ( $^{\circ}\text{C}$ )	29 31
REL.Humid. ( % )	63 67
AC Supply ( Volt )	220 220

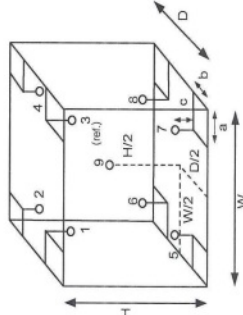

### Probe Installation Details :

a = 10 cm  
b = 10 cm  
c = 10 cm  
D = 0.62 m  
W = 1.2 m  
H = 1.2 m  
Capacity = 0.89 m<sup>3</sup>

Position :	Ref. Std. ID No.:
1	22-18RTD-2/1
2	18RTD-2/2
3	18RTD-2/3
4	18RTD-2/4
5	18RTD-2/5
6	18RTD-2/6
7	18RTD-2/7
8	18RTD-2/8
9 (ref.)	18RTD-2/9

เอกสารไม่ควบคุม

44444444





Cert. No.: 23TM249  
Page : 3 of 3

Function of UUC*	Temperature Source
1. $UUC = 1$	1. $T = 1$
2. $UUC = 2$	2. $T = 2$
3. $UUC = 3$	3. $T = 3$
4. $UUC = 4$	4. $T = 4$
5. $UUC = 5$	5. $T = 5$
6. $UUC = 6$	6. $T = 6$
7. $UUC = 7$	7. $T = 7$
8. $UUC = 8$	8. $T = 8$
9. $UUC = 9$	9. $T = 9$
10. $UUC = 10$	10. $T = 10$
11. $UUC = 11$	11. $T = 11$
12. $UUC = 12$	12. $T = 12$
13. $UUC = 13$	13. $T = 13$
14. $UUC = 14$	14. $T = 14$
15. $UUC = 15$	15. $T = 15$
16. $UUC = 16$	16. $T = 16$
17. $UUC = 17$	17. $T = 17$
18. $UUC = 18$	18. $T = 18$
19. $UUC = 19$	19. $T = 19$
20. $UUC = 20$	20. $T = 20$
21. $UUC = 21$	21. $T = 21$
22. $UUC = 22$	22. $T = 22$
23. $UUC = 23$	23. $T = 23$
24. $UUC = 24$	24. $T = 24$
25. $UUC = 25$	25. $T = 25$
26. $UUC = 26$	26. $T = 26$
27. $UUC = 27$	27. $T = 27$
28. $UUC = 28$	28. $T = 28$
29. $UUC = 29$	29. $T = 29$
30. $UUC = 30$	30. $T = 30$
31. $UUC = 31$	31. $T = 31$
32. $UUC = 32$	32. $T = 32$
33. $UUC = 33$	33. $T = 33$
34. $UUC = 34$	34. $T = 34$
35. $UUC = 35$	35. $T = 35$
36. $UUC = 36$	36. $T = 36$
37. $UUC = 37$	37. $T = 37$
38. $UUC = 38$	38. $T = 38$
39. $UUC = 39$	39. $T = 39$
40. $UUC = 40$	40. $T = 40$
41. $UUC = 41$	41. $T = 41$
42. $UUC = 42$	42. $T = 42$
43. $UUC = 43$	43. $T = 43$
44. $UUC = 44$	44. $T = 44$
45. $UUC = 45$	45. $T = 45$
46. $UUC = 46$	46. $T = 46$
47. $UUC = 47$	47. $T = 47$
48. $UUC = 48$	48. $T = 48$
49. $UUC = 49$	49. $T = 49$
50. $UUC = 50$	50. $T = 50$
51. $UUC = 51$	51. $T = 51$
52. $UUC = 52$	52. $T = 52$
53. $UUC = 53$	53. $T = 53$
54. $UUC = 54$	54. $T = 54$
55. $UUC = 55$	55. $T = 55$
56. $UUC = 56$	56. $T = 56$
57. $UUC = 57$	57. $T = 57$
58. $UUC = 58$	58. $T = 58$
59. $UUC = 59$	59. $T = 59$
60. $UUC = 60$	60. $T = 60$
61. $UUC = 61$	61. $T = 61$
62. $UUC = 62$	62. $T = 62$
63. $UUC = 63$	63. $T = 63$
64. $UUC = 64$	64. $T = 64$
65. $UUC = 65$	65. $T = 65$
66. $UUC = 66$	66. $T = 66$
67. $UUC = 67$	67. $T = 67$
68. $UUC = 68$	68. $T = 68$
69. $UUC = 69$	69. $T = 69$
70. $UUC = 70$	70. $T = 70$
71. $UUC = 71$	71. $T = 71$
72. $UUC = 72$	72. $T = 72$
73. $UUC = 73$	73. $T = 73$
74. $UUC = 74$	74. $T = 74$
75. $UUC = 75$	75. $T = 75$
76. $UUC = 76$	76. $T = 76$
77. $UUC = 77$	77. $T = 77$
78. $UUC = 78$	78. $T = 78$
79. $UUC = 79$	79. $T = 79$
80. $UUC = 80$	80. $T = 80$
81. $UUC = 81$	81. $T = 81$
82. $UUC = 82$	82. $T = 82$
83. $UUC = 83$	83. $T = 83$
84. $UUC = 84$	84. $T = 84$
85. $UUC = 85$	85. $T = 85$
86. $UUC = 86$	86. $T = 86$
87. $UUC = 87$	87. $T = 87$
88. $UUC = 88$	88. $T = 88$
89. $UUC = 89$	89. $T = 89$
90. $UUC = 90$	90. $T = 90$
91. $UUC = 91$	91. $T = 91$

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Uncertainty (± °C)	Coverage Factor <i>k</i>
20.0	20.0	19.3	0.32	0.57	1.0	0.60	2

Calibration Point ( °C )	Measured Temperature ( °C )								
	Position								
	1	2	3	4	5	6	7	8	9 (ref.)
20.0	20.086	19.916	20.386	19.976	19.973	19.838	19.837	19.821	19.949

**Temperature stability :** One-half of the greatest maximum difference of measured temperature at any one sensor.

temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

UUC\* : Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage

-000-

# เอกสารไม่ควบคุม

a 1149512



NSC-TISI-TIS17025  
CALIBRATION 0008

Cert. No.: 23TM192  
Page: 1 of 3

# Certificate of Calibration

Equipment : Incubator

**Manufacturer :**

Model : BD 53 E2

Serial No. : 13-07343

ID No.: UAE.MIC.005/2558

Submitted by :  
United Analyst and Engineering Consultant Co., Ltd.

3 Soi Udomsuk 41, Sukhumvit Road,

Banachak, Phrakhanong.

Bangkok 10260

**Location :** Microbiology Laboratory

Received Order : 15 February 2023

Calibration Date : 15 February 2023

Ambient Temperature :  $(26 \pm 10)^\circ\text{C}$ Relative Humidity :  $(50 \pm 30) \%$ 

Calibrated by : Suwit Imjai

Approved by :

( ) Pornthippa Tameyakul

( ) Malee Butkruea

Issue Date : 24 February 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services

เขียน  
Services.  
เอกสารไม่ควบคุม







TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES  
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250  
TEL. 0-2717-3000-29 FAX. 0-2719-9484



NSC-TS1-TS17025  
CALIBRATION 0008

Cert. No.: 23TM193  
Page : 1 of 3

## Certificate of Calibration

**Equipment :** Water Bath  
**Manufacturer :** Memmert  
**Model :** WNE 14  
**Serial No. :** L416.0606  
**ID No. :** UAE.MIC.002/2560

**Submitted by :** United Analyst and Engineering Consultant Co., Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260

**Location :** Microbiology Laboratory

**Received Order :** 15 February 2023  
**Calibration Date :** 15 February 2023  
**Ambient Temperature :** ( 26 ± 10 ) °C  
**Relative Humidity :** ( 50 ± 30 ) %

**Calibrated by :** Suwit Imjai

**Approved by :**

( ) Pornthippa Tameyakul  
( ) Malee Buikruea

**Issue Date :** 24 February 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services

เอกสารไม่ควบคุม



**Equipment :** Water Bath  
**Condition As-Received :** Used Item  
**Reference :** 2302-02950C-2  
**Procedure Used :-**

Calibration were conducted using in-house calibration procedure CP-OT04 according to direct measurement method with Data Acquisition which connected with Industrial Platinum Resistance Thermometer ( IPRT ).

The temperature scale used was based on ITS-90.

### Condition of this result of calibration

1. Reference standard instrument:-

**Instrument** **Model** **Serial No.** **Cert. No.** **Due Date**  
1 ) Data Acquisition 34972A MY59003411 22LM165 26 Nov 2023


2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

**Result of Calibration :-** ( \* ) Without Adjustment

**Function of UUC\* :** Temperature Source

	Environmental		AC Voltage Supply ( Volt )
	( °C )	( %R.H. )	
Beginning of Calibration	22	65	231
Finished of Calibration	23	61	231

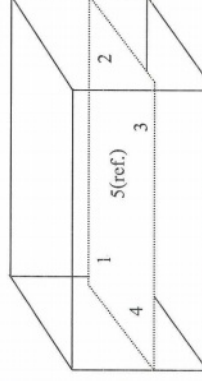

Front

Position :	Ref. Std. ID No.:
1	4804539-001
2	4804539-002
3	4804539-003
4	4804539-004
5(ref.)	4804539-005

Cert. No.: 23TM193  
Page : 2 of 3

เอกสารไม่ควบคุม



Equipment : Water Bath  
Condition As-Received : Used Item  
Reference : 2302-0295OC-2  
Result of Calibration :- ( \* ) Without Adjustment  
Function of UUC\* : Temperature Source

Cert. No.: 23TM193

Page : 3 of 3

Calibration point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Average* Standard Reading (°C)				
			Position				
			1	2	3	4	5 (ref.)
44.5	44.5	44.5	44.453	44.437	44.428	44.477	44.459

Calibration point (°C)	Uniformity (°C)	Stability (± °C)	Uncertainty (± °C)	Coverage Factor k
44.5	0.079	0.038	0.15	2

**Average\*** : The average of 30 values in each position.

**Uniformity** : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

**Stability** : One-half of the greatest maximum difference of measured temperature at any one probe.

**UUC\*** : Unit Under Calibration

**Note** : The reported uncertainty of measurement was included stability and excluded uniformity.

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor  $k$ , providing a level of confidence of approximately 95 %.

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES  
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250  
TEL.: 0-2717-3000-27 FAX: 0-2719-9484



NSC-TS&T-171825  
CALIBRATION 0008

Cert. No.: 22TM1121

Page.: 1 of 3

## Certificate of Calibration

Equipment : Autoclave  
Manufacturer : ALP  
Model : CL-40L  
Serial No. : 807298  
ID No. : UAE.MIC.019/2560  
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road,  
Bangchak, Phrakhanong,  
Bangkok 10260  
301 Room  
Location :  
Received Order : 11 July 2022  
Calibration Date : 11 July 2022  
Ambient Temperature : ( 26 ± 10 ) °C  
Relative Humidity : ( 50 ± 30 ) %  
Calibrated by : Preecha Hlahib

Approved by :

( / ) Pornthippa Tameyakul  
( / ) Malee Butkruea  
( ) Suwit Injai

Issue Date :

18 July 2022

The Uncertainties are for a confidence probability of approximately 95 %

This certificate may not be reproduced other than in full, except with the prior written

Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services

เอกสารไม่ควบคุม

เอกสารไม่ควบคุม





Equipment : Autoclave  
Condition As-Received : Used Item  
Reference : 2207-0245OC-7  
Procedure Used :-

Cert. No.: 22TM1121  
Page.: 2 of 3

Calibration were conducted using in-house calibration procedure CP-OT03 according to direct measurement method with Data Acquisition which connected with Thermocouple Type T  
The temperature scale used was based on ITS-90.

#### Condition of this result of calibration

- Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1 ) Data Acquisition	34970A	MY44060450	22LM46	28 Mar 2023
  - This certificate is valid only to the item calibrated on date and place of calibration.
  - This certification is traceable to the International System of Unit.
  - This result of calibration covers laboratory autoclaves for the sterilization of goods and material which could be infected with organisms categorized as Hazard Group 1, 2 and 3\*\*
- (\*\* = Categorization of pathogens according to hazard and categories of containment, second edition, 1990 )  
It does not cover autoclaves for use with material infect with organisms in Hazard Group 4, for which complete containment and sterilization of infected condensate is considered to be essential.
- This result of calibration does not apply to sterilizers or disinfectors used for medical, dental, pharmaceutical or veterinary purposes which are directly concerned with patient care, or those used for fabrics subjected to sterilization which are required to be dry at the end of cycle.

Result of Calibration :- ( \* ) Without Adjustment  
Function of UUC\* : Temperature Source



Environmental		Ref. Std. ID No.:	
Beginning of Calibration	( °C )	( %R.H. )	( Volt )
	29	49	220
Finished of Calibration	( °C )	( %R.H. )	( Volt )
	32	48	220
Position	Description		
1 =	Center of chamber		
2 =	Temperature sensor		
3 =	Exhaust port		

เอกสารไม่ควบคุม



Equipment : Autoclave  
Condition As-Received : Used Item  
Reference : 2207-0245OC-7  
Result of Calibration :- ( \* ) Without Adjustment

Cert. No.: 22TM1121  
Page.: 3 of 3

Operating parameter Set : Temperature = 115 °C  
Sterilization period = 15 minute

UUC* Setting ( °C )	UUC* Reading ( °C )	Position	Average* Standard Reading ( °C )	Stability ( ± °C )	Pressure Reading ( MPa )	Uncertainty ( ± °C )	Coverage Factor k
116	116	1	116.523	0.14	0.08	0.90	2
		2	116.566				
		3	116.440				

Operating parameter Set : Temperature = 121 °C  
Sterilization period = 30 minute

UUC* Setting ( °C )	UUC* Reading ( °C )	Position	Average* Standard Reading ( °C )	Stability ( ± °C )	Pressure Reading ( MPa )	Uncertainty ( ± °C )	Coverage Factor k
122	122	1	122.503	0.19	0.12	0.91	2
		2	122.637				
		3	122.558				

Average\* : The average of 30 values in each position.

Stability : One-half of the greatest maximum difference of measured temperature at any one probe.

UUC\* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่ควบคุม



## Certificate of Calibration

**Equipment:** Balance  
**Model:** PX623  
**Serial No. (or ID.):** C236754745  
**Manufacturer:** Ohaus  
**Condition:** New

Certificate No.: C01223732

Issued Date: 09 December 2022

Job No.: KSPR2215576

Page: 1 of 2

**Customer:** United Analyst and Engineering Consultant Co., Ltd.  
3 Soi Udomsuk 41, Sukhumvit Road, Bangchak Sub-District,  
Phrakhanong District, Bangkok, THAILAND 10260

**Environment Condition:** Temperature 26 °C ± 0.5 °C  
Humidity 53 %RH ± 3.9 %RH

**Calibration Place:** United Analyst and Engineering Consultant Co., Ltd. (301 Microbiology Room)  
3 Soi Udomsuk 41, Sukhumvit Road, Bangchak Sub-District,  
Phrakhanong District, Bangkok, THAILAND 10260

**Calibration By:** Mr. Adisai Maknoi  
**Calibration Date:** 09 December 2022  
**The Method used:** In-house method, CAL-WI-47, based on UKAS Lab 14  
**Traceability:** This certificate is traceable to the SI Units maintained by National Institute of Metrology (NIMT), Thailand through DKSH Technology Co., Ltd. Certificate No. C02221765

*Ad.*

(Mr. Adisai Maknoi)

Person in charge

This certificate is issued the units of measurement according to the International System of Units (SI). It provides traceability of measurement to International or national standard or other recognized national standard laboratories.  
The measurement uncertainty stated is the expanded uncertainty which is obtained from the standard uncertainty multiplied by the coverage factor ( $k=2$ ) to provide a level of confidence of approximately 95%. It is determined in accordance with the Guide to Expression of Uncertainty in Measurement (GUM).  
These results may be affected by deviations from specified conditions. The results relate only to the items tested, calibrated or sampled. The report shall not be reproduced except in full without approval of DKSH Technology Limited.

บริษัท ดีเคเอส อีเซีย จำกัด  
DKSH Technology Limited  
2533 ถนนสุขุมวิท แขวงคลองเตย เขตคลองเตย กรุงเทพมหานคร 10260  
2533 Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260  
Phone: +66 2639 7000 Email: info.calibration@dksh.com Website: www.dksh.com/scientific-thailand

Delivering Growth – in Asia and Beyond.

CAL-FM-C01-14: 12 Sep 2022



Certificate No.: C01223732

Page: 2 of 2

## Calibration Results: Without Adjustment

Eccentric Error: Weight to be 1/3 or 1/2 of Maximum capacity, taken from the center of the pan as a zero reference.

Nominal Test Value	Reference Points (g)				
	A	B	C	D	E
-	-	0.000	0.000	0.000	0.000

Repeatability: Determination of the standard deviation of weighing balance., Readability 0.001 (g)

Nominal test value (g)	Standard Deviation
50	0.0004
500	0.0005

Error of Indication from nominal or conventional mass value., Readability 0.001 (g)

Nominal Value (g)	Conventional Mass (g)	Displayed Value (g)	Error of Indication (g)	Uncertainty (g)	k
1	1.0000	1.000	0.000	0.0010	2.03
5	5.0001	5.000	0.000	0.0010	2.03
10	10.0001	10.000	0.000	0.0010	2.03
20	20.0001	20.000	0.000	0.0010	2.03
50	50.0001	50.000	0.000	0.0010	2.03
100	100.0001	100.000	0.000	0.0011	2.03
200	200.0004	200.000	0.000	0.0011	2.02
300	300.0005	300.000	-0.001	0.0013	2.01
400	400.0008	400.001	0.000	0.0014	2.01
500	500.0003	500.000	0.000	0.0017	2.00
600	600.0004	600.000	0.000	0.0019	2.00

The End of Certificate

*Rungrod*

(Mr. Rungrod Jenkitrakulchai)

Authorized signatory

เอกสารไม่ควบคุม

Delivering Growth – in Asia and Beyond.

CAL-FM-C01-14: 12 Sep 2022



บริษัท ดีเคเอส อีเซีย จำกัด  
DKSH Technology Limited  
2533 ถนนสุขุมวิท แขวงคลองเตย เขตคลองเตย กรุงเทพมหานคร 10260  
2533 Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260  
Phone: +66 2639 7000 Email: info.calibration@dksh.com Website: www.dksh.com/scientific-thailand

Delivering Growth – in Asia and Beyond.

CAL-FM-C01-14: 12 Sep 2022

เอกสารไม่ควบคุม



## Statements of conformity:

This conformity certificate documents the validity of the following statements of conformity based on the measurement results of corresponding calibration certificate:

The error of indication determined during calibration are under given measurement and environmental conditions and considering the expanded measurement uncertainty (coverage probability 95%) within the specification. The given measurement uncertainty already includes other all effects by according to the standard method, UKAS Lab14. Therefore, those parameters have not been assessed separately.

### Tolerance and Decision rules:

Assessment of the conformity of the measurement device are done based on direct comparison of the relevant measurement results with the tolerances and decision rule are prescribed by the customer.

Decision rule : ☐ Choice A Binary Statement for Simple Acceptance Rule ( $w = 0$ ), Specific Risk < 50% PFA.

☒ Choice B Non-binary statement with guard band ( $w = 1$  U), Pass or Fail Specific Risk < 2.5% PFA and Condition Pass or Condition Fail Specific Risk < 50% PFA.

☐ Choice C Customer defined, Customers may define arbitrary multiple of  $r$  to have applied as guard band ( $w = r$  U).

; PFA – Probability of False Accept



Authorized signatory

The validity of the statements of conformity cannot be guaranteed for different places of use, environmental conditions or improper use.

## The End of Statements of conformity


## Statements of conformity:

### Without Adjustment

Readability: 0.001 g

Nominal Value g	Error of indication g	Guard band (w) g	Tolerance (±) g	Conformity
1	0.000	0.0010	0.002	Pass
5	0.000	0.0010	0.010	Pass
10	0.000	0.0010	0.020	Pass
20	0.000	0.0010	0.040	Pass
50	0.000	0.0010	0.100	Pass
100	0.000	0.0011	0.200	Pass
200	0.000	0.0011	0.400	Pass
300	-0.001	0.0013	0.600	Pass
400	0.000	0.0014	0.800	Pass
500	0.000	0.0017	1.000	Pass
600	0.000	0.0019	1.200	Pass



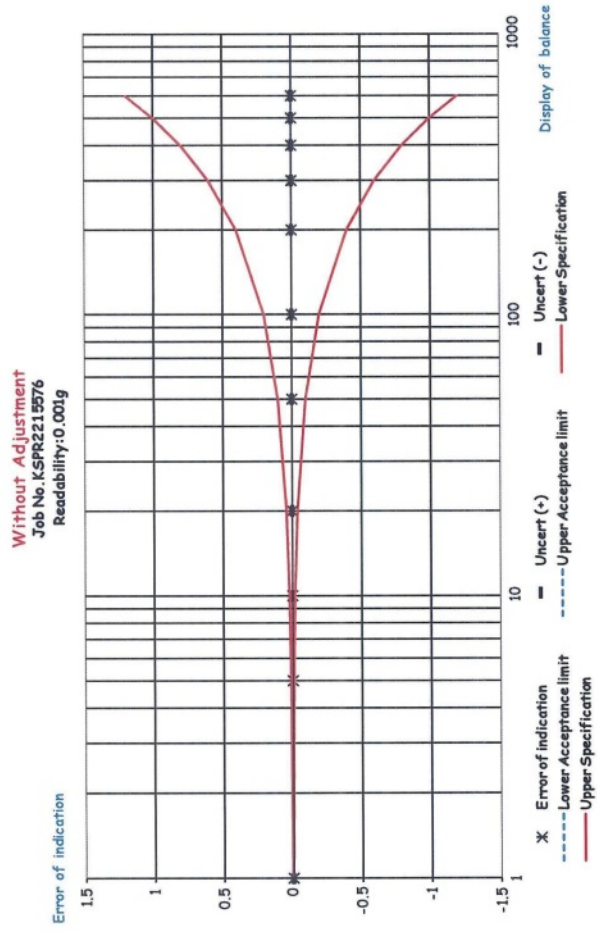

เอกสารไม่ควบคุม

ภาคผนวก ฉ  
หนังสือขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน



- ๒ -

ที่ อก ๐๓๑๐(๑)/ ๖ ๐ ๒ ๘

กรมโรงงานอุตสาหกรรม  
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท  
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๒๒ มีนาคม ๒๕๖๖

เรื่อง เปลี่ยนแปลงบุคลากรและสารมลพิษที่วิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท ยูโนเด็ค แอนนาลิสต์ แอนด์ เอ็นจิเนียริ่ง คอนซัลแตนท์ จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน  
ลงวันที่ ๓๐ มกราคม ๒๕๖๖

สิ่งที่ส่งมาด้วย เอกสารแนบท้ายหนังสือเปลี่ยนแปลงบุคลากรและสารมลพิษที่วิเคราะห์

บริษัท ยูโนเด็ค แอนนาลิสต์ แอนด์ เอ็นจิเนียริ่ง คอนซัลแตนท์ จำกัด จำนวน ๒ แผ่น

ตามหนังสือที่อ้างถึง บริษัท ยูโนเด็ค แอนนาลิสต์ แอนด์ เอ็นจิเนียริ่ง คอนซัลแตนท์ จำกัด  
ห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ๖-๑๔๕ สถานที่ตั้งเลขที่ ๓ ซอยอุดมสุข ๔๑ ถนนสุขุมวิท  
แขวงบางจาก เขตพระโขนง กรุงเทพมหานคร ขอเปลี่ยนแปลงบุคลากรและสารมลพิษที่วิเคราะห์  
ความละเอียดแล้วแล้ว นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว มีความเห็นดังนี้

๑. ให้ยกเลิกเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๒ ราย

๑) [REDACTED] ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๐๑๖

๒) [REDACTED] ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๐๕๗

๒. ให้เพิ่มเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๑๑ ราย

๑) [REDACTED] ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๑๓๒

๒) [REDACTED] ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๑๓๓

๓) [REDACTED] ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๑๓๔

๔) [REDACTED] ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๑๓๕

๕) [REDACTED] ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๑๓๖

๖) [REDACTED] ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๑๓๗

๗) [REDACTED] ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๑๓๘

๘) [REDACTED] ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๑๓๙

๙) [REDACTED] ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๑๔๐

๑๐) [REDACTED] ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๑๔๑

๑๑) [REDACTED] ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๑๔๒

๓. ให้เพิ่มขอบข่ายสารมลพิษ

อนึ่ง...

อนึ่ง หนังสือฉบับนี้จะมีผลใช้บังคับเมื่อพ้นกำหนดสิบวันนับแต่วันพ้นจากอำนาจของเจ้าพนักงาน  
ที่ อก ๐๓๑๐(๑)/๑๘๘๙ ลงวันที่ ๙ กุมภาพันธ์ ๒๕๖๕ คือในวันที่ ๒ กุมภาพันธ์ ๒๕๖๘ ทั้งนี้ สามารถยื่นคำขอ  
ผ่านระบบอิเล็กทรอนิกส์ได้ที่หน้าเว็บไซต์กรมโรงงานอุตสาหกรรม ตาม QR Code พ้ายหนังสือฉบับนี้

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ



ผู้อำนวยการกองวิจัยและเตือนภัยมลพิษโรงงาน  
ปฏิบัติการตามแผนป้องกันมลพิษโรงงานอุตสาหกรรม



ยื่นคำขอผ่านระบบอิเล็กทรอนิกส์

กองวิจัยและเตือนภัยมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓๕ โทรสาร ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๑๙๙

ไปรษณีย์อิเล็กทรอนิกส์ sarabang@dlw.mail.go.th



"อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว"



- ๒ -

เอกสารแนบท้ายหนังสือเปลี่ยนแปลงบุคลากรและสารมลพิษที่วิเคราะห์

บริษัท ยูโนเด็ค แอนนาลิสต์ แอนด์ เอ็นจิเนียริ่ง คอนซัลแตนท์ จำกัด เลขทะเบียน ๖-๑๔๕

ที่ อก ๐๓๑๐(๑)/ ๖ ๐ ๒ ๘

ลงวันที่ ๒๒ มีนาคม ๒๕๖๖

ขอบข่ายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๑๖ รายการ

เดิม จำนวน ๑๖ รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Benzene	Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method <sup>(1,2)</sup>
2	Carbon tetrachloride	Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method <sup>(1,2)</sup>
3	1,2-Dichloroethane	Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method <sup>(1,2)</sup>
4	1,1-Dichloroethylene	Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method <sup>(1,2)</sup>
5	cis-1,2-Dichloroethylene	Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method <sup>(1,2)</sup>
6	trans-1,2-Dichloroethylene	Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method <sup>(1,2)</sup>
7	Ethylbenzene	Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method <sup>(1,2)</sup>
8	Methylene chloride	Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method <sup>(1,2)</sup>
9	Styrene	Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method <sup>(1,2)</sup>
10	Tetrachloroethylene	Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method <sup>(1,2)</sup>
11	Toluene	Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method <sup>(1,2)</sup>
12	Trichloroethylene	Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method <sup>(1,2)</sup>
13	m-Xylene	Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method <sup>(1,2)</sup>
14	o-Xylene	Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method <sup>(1,2)</sup>
15	p-Xylene	Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method <sup>(1,2)</sup>
16	Xylene (Total)	Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method <sup>(1,2)</sup>

เอกสารอ้างอิง...

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ กองวิจัยและเตือนภัยมลพิษโรงงาน กรมโรงงานอุตสาหกรรม โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓๕





ที่ อก ๐๓๑๐(๑)/ ๑๕๕๕ ๓

กรมโรงงานอุตสาหกรรม  
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท  
เขตราชเทวี กรุงเทพฯ ๑๐๕๐๐

๒๕ ตุลาคม ๒๕๖๕

เรื่อง เปลี่ยนแปลงบุคลากรของห้องปฏิบัติการวิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท ยูโนเด็ค แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารเคมีของห้องปฏิบัติการวิเคราะห์เอกชน  
ลงวันที่ ๑๘ ตุลาคม ๒๕๖๕

ตามหนังสือที่อ้างถึง บริษัท ยูโนเด็ค แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด  
ห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ๖-๑๕๕๕ สถานที่ตั้งเลขที่ ๓ ซอยอุดมสุข ๔๑ ถนนสุขุมวิท แขวงบางจาก  
เขตพระโขนง กรุงเทพมหานคร ขอเปลี่ยนแปลงบุคลากรของห้องปฏิบัติการวิเคราะห์ ความละเอียดแจ้งแล้ว นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว มีความเห็นดังนี้

๑. ให้อยกเลิกเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๔ ราย

- |               |                             |
|---------------|-----------------------------|
| ๑) [REDACTED] | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๐๐๒ |
| ๒) [REDACTED] | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๐๒๕ |
| ๓) [REDACTED] | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๐๔๕ |
| ๔) [REDACTED] | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๐๔๗ |

๒. ให้เพิ่มผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๒ ราย

- |               |                             |
|---------------|-----------------------------|
| ๑) [REDACTED] | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๐๔๑ |
| ๒) [REDACTED] | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๐๔๒ |

๓. ให้เพิ่มเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๑๒ ราย

- |                |                             |
|----------------|-----------------------------|
| ๑) [REDACTED]  | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๑๒๐ |
| ๒) [REDACTED]  | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๑๒๑ |
| ๓) [REDACTED]  | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๑๒๒ |
| ๔) [REDACTED]  | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๑๒๓ |
| ๕) [REDACTED]  | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๑๒๔ |
| ๖) [REDACTED]  | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๑๒๕ |
| ๗) [REDACTED]  | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๑๒๖ |
| ๘) [REDACTED]  | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๑๒๗ |
| ๙) [REDACTED]  | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๑๒๘ |
| ๑๐) [REDACTED] | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๑๒๙ |
| ๑๑) [REDACTED] | [REDACTED]                  |
| ๑๒) [REDACTED] | [REDACTED]                  |

อนึ่ง หนังสือฉบับนี้...



ที่ อก ๐๓๑๐(๑)/ ๑๒ ๑๕๕ ๓

กรมโรงงานอุตสาหกรรม  
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท  
เขตราชเทวี กรุงเทพฯ ๑๐๕๐๐

๐๑ กันยายน ๒๕๖๕

เรื่อง เปลี่ยนแปลงบุคลากรของห้องปฏิบัติการวิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท ยูโนเด็ค แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารเคมีของห้องปฏิบัติการวิเคราะห์เอกชน  
ลงวันที่ ๒๖ สิงหาคม ๒๕๖๕

ตามหนังสือที่อ้างถึง บริษัท ยูโนเด็ค แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด  
ห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ๖-๑๕๕๕ สถานที่ตั้งเลขที่ ๓ ซอยอุดมสุข ๔๑ ถนนสุขุมวิท แขวงบางจาก  
เขตพระโขนง กรุงเทพมหานคร ขอเปลี่ยนแปลงบุคลากรของห้องปฏิบัติการวิเคราะห์ ความละเอียดแจ้งแล้ว นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว มีความเห็นดังนี้

๑. ให้อยกเลิกเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๘ ราย

- |               |                             |
|---------------|-----------------------------|
| ๑) [REDACTED] | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๐๓๓ |
| ๒) [REDACTED] | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๐๓๕ |
| ๓) [REDACTED] | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๐๔๑ |
| ๔) [REDACTED] | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๐๕๐ |
| ๕) [REDACTED] | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๐๖๓ |
| ๖) [REDACTED] | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๐๘๐ |
| ๗) [REDACTED] | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๐๘๑ |
| ๘) [REDACTED] | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๑๐๘ |

๒. ให้เพิ่มเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๑๑ ราย

- |                |                             |
|----------------|-----------------------------|
| ๑) [REDACTED]  | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๑๐๙ |
| ๒) [REDACTED]  | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๑๑๐ |
| ๓) [REDACTED]  | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๑๑๑ |
| ๔) [REDACTED]  | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๑๑๒ |
| ๕) [REDACTED]  | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๑๑๓ |
| ๖) [REDACTED]  | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๑๑๔ |
| ๗) [REDACTED]  | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๑๑๕ |
| ๘) [REDACTED]  | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๑๑๖ |
| ๙) [REDACTED]  | ทะเบียนเลขที่ ๖-๑๕๕๕-๖-๐๑๑๗ |
| ๑๐) [REDACTED] | [REDACTED]                  |
| ๑๑) [REDACTED] | [REDACTED]                  |

อนึ่ง หนังสือฉบับนี้...

- ๒ -

อนึ่ง หนังสือฉบับนี้จะมีผลใช้บังคับเมื่อได้รับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน  
ที่ อก ๐๓๑๐(๑)/๑๕๕๕ ลงวันที่ ๔ กุมภาพันธ์ ๒๕๖๕ คือในวันที่ ๒ กุมภาพันธ์ ๒๕๖๕ ทั้งนี้ สามารถยื่นคำขอ  
ผ่านระบบอิเล็กทรอนิกส์ได้ที่หน้าเว็บไซต์กรมโรงงานอุตสาหกรรม ตาม QR Code ที่แนบหนังสือฉบับนี้

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ



ผู้อำนวยการกองวิจัยและเตือนภัยมลพิษโรงงาน  
ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม



ยื่นคำขอผ่านระบบอิเล็กทรอนิกส์

กองวิจัยและเตือนภัยมลพิษโรงงาน  
กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ  
โทร. ๐ ๒๕๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕  
โทรสาร ๐ ๒๕๓๐ ๖๓๑๒ ต่อ ๒๑๑๗  
ไปรษณีย์อิเล็กทรอนิกส์ saraban@dlw.mail.go.th



"อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว"



- ๒ -

อนึ่ง หนังสือฉบับนี้จะมีผลใช้บังคับเมื่อได้รับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน  
ที่ อก ๐๓๑๐(๑)/๑๕๕๕ ลงวันที่ ๔ กุมภาพันธ์ ๒๕๖๕ คือในวันที่ ๒ กุมภาพันธ์ ๒๕๖๕ ทั้งนี้ สามารถยื่นคำขอ  
ผ่านระบบอิเล็กทรอนิกส์ได้ที่หน้าเว็บไซต์กรมโรงงานอุตสาหกรรม ตาม QR Code ที่แนบหนังสือฉบับนี้

จึงเรียนมาเพื่อทราบ



ผู้อำนวยการกองวิจัยและเตือนภัยมลพิษโรงงาน  
ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม



ยื่นคำขอผ่านระบบอิเล็กทรอนิกส์

กองวิจัยและเตือนภัยมลพิษโรงงาน  
กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ  
โทร. ๐ ๒๕๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕  
โทรสาร ๐ ๒๕๓๐ ๖๓๑๒ ต่อ ๒๑๑๗  
ไปรษณีย์อิเล็กทรอนิกส์ saraban@dlw.mail.go.th



"อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว"



ที่ อก ๐๓๑๐(๑)/ ๔๗๘ ๗



กรมโรงงานอุตสาหกรรม  
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท  
เขตราชเทวี กรุงเทพฯ ๑๐๕๐๐

๒๑ เมษายน ๒๕๖๕

เรื่อง เปลี่ยนแปลงบุคลากรของห้องปฏิบัติการวิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท ยูโนเด็ค แอนนาลิสต์ แอนด์ เอ็นจิเนียริ่ง คอนซัลแตนท์ จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารเคมีของห้องปฏิบัติการวิเคราะห์ ความละเอียดแจ้งแล้ว นั้น  
ลงวันที่ ๓๐ มีนาคม ๒๕๖๕

ตามหนังสือที่อ้างถึง บริษัท ยูโนเด็ค แอนนาลิสต์ แอนด์ เอ็นจิเนียริ่ง คอนซัลแตนท์ จำกัด  
ห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ๗-๑๕๕ สถานที่ตั้งเลขที่ ๓ ซอยอุดมสุข ถนนสุขุมวิท แขวงบางจาก  
เขตพระโขนง กรุงเทพมหานคร ขอเปลี่ยนแปลงบุคลากรของห้องปฏิบัติการวิเคราะห์ ความละเอียดแจ้งแล้ว นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว มีความเห็นดังนี้

๑. ให้ยกเลิกควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๒ ราย

๑) [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๐๕

๒) [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๐๖

๒. ให้ยกเลิกเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๒ ราย

๑) [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-จ-๐๐๖๔

๒) [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-จ-๐๐๘๔

๓. ให้เพิ่มเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๒ ราย

๑) [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-จ-๐๑๐๗

๒) [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-จ-๐๑๐๘

อนึ่ง หนังสือฉบับนี้จะตามต่ออายุพร้อมหนังสือขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน  
ที่ อก ๐๓๑๐(๑)/๑๘๗๙ ลงวันที่ ๙ กุมภาพันธ์ ๒๕๖๕ คือในวันที่ ๒ กุมภาพันธ์ ๒๕๖๘ ทั้งนี้ สามารถยื่นคำขอ  
ผ่านระบบอิเล็กทรอนิกส์ได้ที่หน้าเว็บไซต์กรมโรงงานอุตสาหกรรม ตาม QR Code หายหนังสือฉบับนี้

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

(นางจุฬารัตน์ เจริญกิจ  
ผู้อำนวยการกองทะเบียนและใบอนุญาตประกอบกิจการโรงงาน  
ปฏิบัติการและควบคุมโรงงานอุตสาหกรรม)



ยื่นคำขอผ่านระบบอิเล็กทรอนิกส์

กองวิจัยและเตือนภัยมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและ

โทร. ๐ ๒๕๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕ โทรสาร

ไปรษณีย์อิเล็กทรอนิกส์ saraban@dlw.mail.go.th



"อุตสาหกรรมก้าวไกล ประเทศไทยก้าวร้าว ร่วมกันพัฒนา อุตสาหกรรมสีเขียว"



ที่ อก ๐๓๑๐(๑)/ ๑๘๗ ๘

กรมโรงงานอุตสาหกรรม  
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท  
เขตราชเทวี กรุงเทพฯ ๑๐๕๐๐

๐๙ กุมภาพันธ์ ๒๕๖๕

เรื่อง ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

เรียน กรรมการผู้จัดการ บริษัท ยูโนเด็ค แอนนาลิสต์ แอนด์ เอ็นจิเนียริ่ง คอนซัลแตนท์ จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารเคมีของห้องปฏิบัติการวิเคราะห์เอกชน  
ลงวันที่ ๒๗ ธันวาคม ๒๕๖๔

สิ่งที่ส่งมาด้วย ๑. รายชื่อผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๔๐ ราย

๒. รายชื่อเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๑๐๖ ราย

๓. ขอบข่ายสารเคมีที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม

ตามหนังสือที่อ้างถึง บริษัท ยูโนเด็ค แอนนาลิสต์ แอนด์ เอ็นจิเนียริ่ง คอนซัลแตนท์ จำกัด  
ขอต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ๗-๑๕๕ สถานที่ตั้งเลขที่ ๓  
ซอยอุดมสุข ๔๑ ถนนสุขุมวิท แขวงบางจาก เขตพระโขนง กรุงเทพมหานคร ต่อกรมโรงงานอุตสาหกรรม นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้บริษัท ยูโนเด็ค แอนนาลิสต์ แอนด์ เอ็นจิเนียริ่ง  
คอนซัลแตนท์ จำกัด ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน โดยมีองค์ประกอบดังนี้

ก. ผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๔๐ ราย ตามสิ่งที่ส่งมาด้วย ๑

ข. เจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๑๐๖ ราย ตามสิ่งที่ส่งมาด้วย ๒

ค. ขอบข่ายสารเคมีที่ได้รับขึ้นทะเบียนไว้วิเคราะห์ในน้ำเสีย น้ำใต้ดิน อากาศเสีย สิ่งปฏิกูล  
หรือวัสดุที่ไม่ใช้แล้ว และดิน ตามสิ่งที่ส่งมาด้วย ๓

หนังสือฉบับนี้จะตามต่ออายุในวันที่ ๒ กุมภาพันธ์ ๒๕๖๘ หากประสงค์จะต่ออายุหนังสือ  
รับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ให้ยื่นคำขอต่ออายุพร้อมเอกสารประกอบคำขอต่อ  
กรมโรงงานอุตสาหกรรมภายใน ๓๐ วัน ก่อนวันสิ้นสุดของหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน  
ทั้งนี้ สามารถยื่นคำขอผ่านระบบอิเล็กทรอนิกส์ได้ที่หน้าเว็บไซต์กรมโรงงานอุตสาหกรรม ตาม QR Code หาย  
หนังสือฉบับนี้

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

(นางจุฬารัตน์ เจริญกิจ  
ผู้อำนวยการกองทะเบียนและใบอนุญาตประกอบกิจการโรงงาน  
ปฏิบัติการและควบคุมโรงงานอุตสาหกรรม)



ยื่นคำขอผ่านระบบอิเล็กทรอนิกส์

กองวิจัยและเตือนภัยมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและ

โทร. ๐ ๒๕๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕ โทรสาร

ไปรษณีย์อิเล็กทรอนิกส์ saraban@dlw.mail.go.th

ไปรษณีย์อิเล็กทรอนิกส์ saraban@dlw.mail.go.th

สิ่งที่ส่งมาด้วย ๑

เอกสารแนบท้ายหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท ยูโนเด็ค แอนนาลิสต์ แอนด์ เอ็นจิเนียริ่ง คอนซัลแตนท์ จำกัด เลขทะเบียน ๗-๑๕๕

ที่ อก ๐๓๑๐(๑)/ ๑๘๗ ๘

ลงวันที่ ๐๙ กุมภาพันธ์ ๒๕๖๕

ก. ผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๔๐ ราย

๑) [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๐๑

๒) [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๐๒

๓) [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๐๓

๔) [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๐๔

๕) [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๐๕

๖) [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๐๖

๗) [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๐๗

๘) [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๐๘

๙) [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๐๙

๑๐ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๑๐

๑๑ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๑๑

๑๒ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๑๒

๑๓ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๑๓

๑๔ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๑๔

๑๕ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๑๕

๑๖ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๑๖

๑๗ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๑๗

๑๘ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๑๘

๑๙ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๑๙

๒๐ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๒๐

๒๑ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๒๑

๒๒ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๒๒

๒๓ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๒๓

๒๔ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๒๔

๒๕ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๒๕

๒๖ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๒๖

๒๗ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๒๗

๒๘ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๒๘

๒๙ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๒๙

๓๐ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๓๐

๓๑ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๓๑

๓๒ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๓๒

๓๓ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๓๓

๓๔ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๓๔

๓๕ [REDACTED] ทะเบียนเลขที่ ๗-๑๕๕-ค-๐๐๓๕

ผู้อำนวยการกองวิจัยและเตือนภัยมลพิษโรงงาน  
ปฏิบัติการและควบคุมโรงงานอุตสาหกรรม

๓๖) นายศุภณัฐ...

எ) ௭௦)  
எ) ௭௧)  
எ) ௭௨)  
எ) ௭௩)  
எ) ௭௪)

ทะเบียนเลขที่ 7-๓๔๕-ค-๐๐๓๖  
ทะเบียนเลขที่ 7-๓๔๕-ค-๐๐๓๗  
ทะเบียนเลขที่ 7-๓๔๕-ค-๐๐๓๘  
ทะเบียนเลขที่ 7-๓๔๕-ค-๐๐๓๙  
ทะเบียนเลขที่ 7-๓๔๕-ค-๐๐๔๐

ผู้ชำนาญการกองวิจัยและพัฒนาเทคโนโลยีโรงงาน  
ปฏิบัติการราชการแบบอินทิเกรตโรงงานอุตสาหกรรม

UNITED ANALYST AND ENGINEERING  
CONSULTANT COMPANY LIMITED

## สามเณรอุกตอง

เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกสาร  
บริษัท ยูนิเค็ด แอนนาลิสต์ แอนด์ เอ็นจิเนียริ่ง คอนซัลแตนท์ จำกัด เลขทะเบียน ว-๓๔๔  
ที่ อก ๐๓๐๑(๑)/ ๑๔๗/ ๙ ลงวันที่ ๐๙ กุมภาพันธ์ ๒๕๖๕

ข. เจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๑๐๖ ราย

၁) :  
 ၂) :  
 ၃) :  
 ၄) :  
 ၅) :  
 ၆) :  
 ၇) :  
 ၈) :  
 ၉) :  
 ၁၀) :  
 ၁၁) :  
 ၁၂) :  
 ၁၃) :  
 ၁၄) :  
 ၁၅) :  
 ၁၆) :  
 ၁၇) :  
 ၁၈) :  
 ၁၉) :  
 ၂၀) :  
 ၂၁) :  
 ၂၂) :  
 ၂၃) :  
 ၂၄) :  
 ၂၅) :  
 ၂၆) :  
 ၂၇) :  
 ၂၈) :  
 ၂၉) :  
 ၃၀) :  
 ၃၁) :  
 ၃၂) :  
 ၃၃) :  
 ၃၄) :  
 ၃၅) :  
 ၃၆) :  
 ၃၇) :  
 ၃၈) :  
 ၃၉) :  
 ၄၀) :  
 ၄၁) :  
 ၄၂) :  
 ၄၃) :  
 ၄၄) :  
 ၄၅) :  
 ၄၆) :  
 ၄၇) :  
 ၄၈) :  
 ၄၉) :  
 ၅၀) :  
 ၅၁) :  
 ၅၂) :  
 ၅၃) :  
 ၅၄) :  
 ၅၅) :  
 ၅၆) :  
 ၅၇) :  
 ၅၈) :  
 ၅၉) :  
 ၆၀) :  
 ၆၁) :  
 ၆၂) :  
 ၆၃) :  
 ၆၄) :  
 ၆၅) :  
 ၆၆) :  
 ၆၇) :  
 ၆၈) :  
 ၆၉) :  
 ၇၀) :  
 ၇၁) :  
 ၇၂) :  
 ၇၃) :  
 ၇၄) :  
 ၇၅) :  
 ၇၆) :  
 ၇၇) :  
 ၇၈) :  
 ၇၉) :  
 ၈၀) :  
 ၈၁) :  
 ၈၂) :  
 ၈၃) :  
 ၈၄) :  
 ၈၅) :  
 ၈၆) :  
 ၈၇) :  
 ၈၈) :  
 ၈၉) :  
 ၉၀) :  
 ၉၁) :  
 ၉၂) :  
 ၉၃) :  
 ၉၄) :  
 ၉၅) :  
 ၉၆) :  
 ၉၇) :  
 ၉၈) :  
 ၉၉) :  
 ၁၀၀) :

๖๓  
 ๖๔  
 ๖๕  
 ๖๖  
 ๖๗  
 ๖๘  
 ๖๙  
 ๗๐  
 ๗๑  
 ๗๒  
 ๗๓  
 ๗๔  
 ๗๕  
 ๗๖  
 ๗๗  
 ๗๘  
 ๗๙  
 ๘๐  
 ๘๑  
 ๘๒  
 ๘๓  
 ๘๔  
 ๘๕  
 ๘๖  
 ๘๗  
 ๘๘  
 ๘๙  
 ๙๐  
 ๙๑  
 ๙๒  
 ๙๓  
 ๙๔  
 ๙๕  
 ๙๖  
 ๙๗  
 ๙๘  
 ๙๙  
 ๑๐๐

ANALYST AND ENGINEERING  
STANT COMPANY LIMITED

ผู้ชี้แจง: นายอรรถพร ฤกษ์พานิชย์ กรรมการผู้จัดการ บริษัท อีสป จำกัด  
บริษัท ปรึกษาการนันทนาการและสันทนาการ จำกัด

[illegible][illegible]

ผู้ชำนาญการกองวิจัยและเฝ้าระวังมลพิษทางอากาศ  
ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม

๗๓) นายอิทธิพงษ์

[illegible][illegible]

ผู้ชำนาญการกองวิจัยและพัฒนาสิ่งแวดล้อมพิษวิทยา  
ปฏิบัติการราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม



เอกสารแนบท้ายหนังสือรับข้ออาชญาบัตรขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน  
บริษัท ยูนิเทค แอนาไลติกส์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด เลขทะเบียน ๖-๑๔๕  
ที่ กอ ๐๓๑๐(๑)/ ๑๘๗๕ ลงวันที่ ๑๔ กุมภาพันธ์ ๒๕๖๕

ขอข่ายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๓๘๗ รายการ

นำสืบ จำนวน 46 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Aldrin	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
2	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>
3	Barium	Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>
4	α-BHC	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
5	β-BHC	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
6	δ-BHC	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
7	γ-BHC	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
8	Biochemical Oxygen Demand	1) 5-Day BOD Test, Azide Modification Method <sup>(4)</sup> 2) 5-Day BOD Test, Membrane Electrode Method <sup>(4)</sup>
9	Cadmium	1) Digestion, Direct Air-Acetylene Flame Method <sup>(4)</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>(4)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>
10	Chemical Oxygen Demand	1) Closed Reflux, Titrimetric Method <sup>(4)</sup> 2) Closed Reflux, Colorimetric Method <sup>(4)</sup> 3) Open Reflux, Titrimetric Method <sup>(4)</sup>
11	Chlordane	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
12	Chromium	1) Digestion, Direct Air-Acetylene Flame Method <sup>(4)</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>(4)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>
13	Color	ADMI Weighted-Ordinate Spectrophotometric Method <sup>(4)</sup>
14	Copper	1) Digestion, Direct Air-Acetylene Flame Method <sup>(4)</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>(4)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>
15	Cyanide	1) Digestion, Direct Air-Acetylene Flame Method <sup>(4)</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>(4)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
16	o,p'-DDT	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
17	4,4'-DDD	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
18	4,4'-DDE	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
19	4,4'-DDT	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
20	Dieldrin	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
21	Endosulfan I	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
22	Endosulfan II	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
23	Endosulfan sulfate	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
24	Endrin	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
25	Endrin aldehyde	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
26	Formaldehyde	Distillation, Colorimetric Method <sup>(3)</sup>
27	Free Chlorine	1) Iodometric Method <sup>(4)</sup> 2) DPD Ferrous Titrimetric Method <sup>(4)</sup>
28	Heptachlor	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
29	Heptachlor Epoxide	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
30	Hexavalent Chromium	1) Colorimetric Method <sup>(4)</sup> 2) Extraction, Direct Air-Acetylene Flame Method <sup>(4)</sup>
31	Lead	1) Digestion, Direct Air-Acetylene Flame Method <sup>(4)</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>(4)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>
32	Manganese	1) Digestion, Direct Air-Acetylene Flame Method <sup>(4)</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>(4)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>
33	Mercury	Digestion, Cold-Vapor Atomic Absorption Spectrometric Method <sup>(4)</sup>
34	Methoxychlor	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
35	Nickel	1) Digestion, Direct Air-Acetylene Flame Method <sup>(4)</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>(4)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>

16 o,p'-DDT...

36 Oil & Grease...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
36	Oil & Grease	1) Liquid-Liquid, Partition-Gravimetric Method <sup>(4)</sup> 2) Soxhlet Extraction Method <sup>(4)</sup>
37	pH	Electrometric Method <sup>(4)</sup>
38	Phenols	1) Distillation, Chloroform Extraction Method <sup>(4)</sup> 2) Distillation, Direct Photometric Method <sup>(4)</sup>
39	Selenium	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>
40	Sulfide	1) Iodometric Method <sup>(4)</sup> 2) Methylene Blue Method <sup>(4)</sup>
41	Temperature	Laboratory and Field Methods <sup>(4)</sup>
42	Total Dissolved Solids	Dried at 180 °C <sup>(4)</sup>
43	Total Kjeldahl Nitrogen	Semi-Micro-Kjeldahl Method <sup>(4)</sup>
44	Total Suspended Solids	Dried at 103-105 °C <sup>(4)</sup>
45	Trivalent Chromium	1) Digestion, Direct Air-Acetylene Flame Method; Colorimetric Method; Calculation <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma Method; Colorimetric Method; Calculation <sup>(4)</sup>
46	Zinc	1) Digestion, Direct Air-Acetylene Flame Method <sup>(4)</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>(4)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>

นำสืบ จำนวน 126 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Acenaphthene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
2	Acetone	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
3	Aldrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
4	Anthracene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
5	Antimony	Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>
6	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>
7	Atrazine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
8	Barium	1) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>
9	Benz(a)anthracene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
10	Benzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
11	Benzo(b)fluoranthene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
12	Benzo(k)fluoranthene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
13	Benzoic acid	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
14	Benzo(a)pyrene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>

4 Anthracene...

15 Benzo(g,h,i)perylene...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
15	Benzo(g,h,i)perylene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
16	Beryllium	Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>
17	Bis(2-chloroethyl)ether	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
18	Bis(2-ethylhexyl)phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
19	Bromodichloromethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
20	Bromoform	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
21	Butanol	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
22	Butyl benzyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
23	Cadmium	1) Digestion, Direct Air-Acetylene Flame Method <sup>(4)</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>(4)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>
24	Carbazole	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
25	Carbon disulfide	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
26	Carbon tetrachloride	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
27	Chlordane	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
28	p-Chloroaniline	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
29	Chlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>

30 Chlorodibromomethane...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
30	Chlorodibromomethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
31	Chloroform	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
32	2-Chlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
33	Chromium	1) Digestion, Direct Air-Acetylene Flame Method <sup>(4)</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>(4)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>
34	Chromium (III)	1) Digestion, Direct Air-Acetylene Flame Method; Colorimetric Method; Calculation <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma Method; Colorimetric Method; Calculation <sup>(4)</sup>
35	Chromium (VI)	1) Colorimetric Method <sup>(4)</sup> 2) Extraction, Air-Acetylene Flame Method <sup>(4)</sup>
36	Chrysene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
37	Cyanide	Distillation, Colorimetric Method <sup>(4)</sup>
38	2,4-D	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
39	DDD	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
40	DDE	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
41	DDT	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>

42 Dibenzo(a,h)anthracene...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
42	Dibenzo(a,h)anthracene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
43	Di-n-butyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
44	1,2-Dichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
45	1,3-Dichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
46	1,4-Dichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
47	3,3'-Dichlorobenzidine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
48	1,1-Dichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
49	1,2-Dichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
50	1,1-Dichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
51	cis-1,2-Dichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
52	trans-1,2-Dichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
53	2,4-Dichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
54	1,2-Dichloropropane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
55	1,3-Dichloropropane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
56	1,3-Dichloropropene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
57	Dieldrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>

58 Diethyl phthalate...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
58	Diethyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
59	2,4-Dimethylphenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
60	2,4-Dinitrophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
61	2,4-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
62	2,6-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
63	Di-n-Octyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
64	Endosulfan	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
65	Endrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
66	Ethylbenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
67	Fluoranthene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
68	Fluorene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
69	Heptachlor	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>

70 Heptachlor epoxide...



ลำดับ	สารมลพิษ	วิธีวิเคราะห์
70	Heptachlor epoxide	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
71	Hexachlorobenzene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
72	Hexachloro-1,3-butadiene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
73	n-Hexane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
74	$\alpha$ -HCH	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
75	$\beta$ -HCH	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
76	$\gamma$ -HCH	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
77	Hexachlorocyclopentadiene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
78	Hexachloroethane	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
79	Indeno(1,2,3-cd)pyrene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
80	Isophorone	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
81	Lead	1) Digestion, Direct Air-Acetylene Flame Method <sup>(4)</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>(4)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>

82 Manganese...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
82	Manganese	1) Digestion, Direct Air-Acetylene Flame Method <sup>(4)</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>(4)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>
83	Mercury	Digestion, Cold-Vapor Atomic Absorption Spectrometric Method <sup>(4)</sup>
84	Methanol	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
85	Methoxychlor	Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup>
86	Methyl bromide	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
87	Methylene chloride	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
88	2-Methylphenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
89	2-Methylnaphthalene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
90	Methyl tert-butyl ether	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
91	Naphthalene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
92	Nickel	1) Digestion, Direct Air-Acetylene Flame Method <sup>(4)</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>(4)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>
93	Nitrobenzene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
94	N-Nitrosodiphenylamine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
95	N-Nitrosodi-n-propylamine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>

96 Polychlorinated Biphenyls...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
96	Polychlorinated Biphenyls - PCB 1016 - PCB 1221 - PCB 1232 - PCB 1242 - PCB 1248 - PCB 1254 - PCB 1260	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
97	Pentachlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
98	pH	Electrometric Method <sup>(4)</sup>
99	Phenanthrene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
100	Phenol	1) Distillation, Chloroform Extraction Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
101	Pyrene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
102	Selenium	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>
103	Silver	Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>
104	Styrene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
105	1,1,2,2-Tetrachloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
106	Tetrachloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
107	Toluene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>

108 Toxaphene...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
108	Toxaphene	1) Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(4)</sup> 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
109	TPH (C <sub>5</sub> - C <sub>9</sub> )	1) Purge and Trap, Gas Chromatographic Method <sup>(11,21)</sup> 2) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(11,25)</sup>
110	TPH (C <sub>10</sub> - C <sub>16</sub> )	Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(9,21)</sup>
111	TPH (C <sub>17</sub> - C <sub>35</sub> )	Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(9,21)</sup>
112	1,2,4-Trichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
113	1,1,1-Trichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
114	1,1,2-Trichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
115	Trichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
116	2,4,5-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
117	2,4,6-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
118	1,3,5-Trimethylbenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
119	Vanadium	Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>
120	Vinyl acetate	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
121	Vinyl chloride	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
122	m-Xylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
123	o-Xylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>

124 p-Xylene...



ลำดับ	สารมลพิษ	วิธีวิเคราะห์
124	p-Xylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
125	Xylene (Total)	Purge and Trap Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
126	Zinc	1) Digestion, Direct Air-Acetylene Flame Method <sup>(4)</sup> 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method <sup>(4)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup>

## อากาศเสีย (ปล่องระบาย) จำนวน 25 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Antimony	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(5)</sup>
2	Arsenic	1) Isokinetic Sampling, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>(3)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(5)</sup>
3	Cadmium	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method <sup>(5)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(5)</sup>
4	Carbon Monoxide	Instrumental Analyzer Method <sup>(5)</sup>
5	Chlorine	Isokinetic Sampling, Ion Chromatographic Method <sup>(5)</sup>
6	Chromium	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method <sup>(5)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(5)</sup>
7	Cobalt	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(5)</sup>
8	Copper	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method <sup>(5)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(5)</sup>
9	Cresol	Absorption Sampling, Gas Chromatographic Method <sup>(4)</sup>

10 Dioxins/Furans...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
10	Dioxins/Furans	Isokinetic Sampling <sup>(4)</sup>
11	Hydrogen Chloride	Isokinetic Sampling, Ion Chromatographic Method <sup>(5)</sup>
12	Hydrogen Fluoride	Isokinetic Sampling, Ion Chromatographic Method <sup>(5)</sup>
13	Hydrogen Sulfide	Absorption Sampling, Iodometric Method <sup>(5)</sup>
14	Lead	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method <sup>(5)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(5)</sup>
15	Manganese	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method <sup>(5)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(5)</sup>
16	Mercury	Isokinetic Sampling, Digestion, Cold-Vapor Atomic Absorption Spectrometric Method <sup>(3)</sup>
17	Nickel	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method <sup>(5)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(5)</sup>
18	Opacity	Ringelmann's Method <sup>(1)</sup>
19	Oxides of Nitrogen	1) Absorption Sampling, Phenoldisulfonic acid Method <sup>(9)</sup> 2) Instrumental Analyzer Method <sup>(5)</sup>
20	Selenium	1) Isokinetic Sampling, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>(3)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(5)</sup>
21	Sulfur Dioxide	1) Absorption Sampling, Barium-Thorin Titrimetric Method <sup>(5)</sup> 2) Instrumental Analyzer Method <sup>(5)</sup>
22	Sulfuric Acid	Isokinetic Sampling, Barium-Thorin Titrimetric Method <sup>(5)</sup>
23	Total Suspended Particulate	Isokinetic Sampling, Gravimetric Method <sup>(5)</sup>
24	Vanadium	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(5)</sup>
25	Xylene	1) Absorption Sampling, Gas Chromatographic Method <sup>(4)</sup> 2) Absorption Sampling, Gas Chromatographic Method <sup>(4)</sup>

สิ่งปฏิกูล...

## สิ่งปฏิกูลหรือวัสดุที่ไม่ใช้แล้ว จำนวน 35 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Aldrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(2,9,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup>
2	Antimony	Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
3	Arsenic	1) Waste Extraction, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>(2,6,15)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(2,6,13)</sup> 3) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>(7,13)</sup> 4) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
4	Barium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(2,6,13)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
5	Beryllium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(2,6,13)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
6	Cadmium	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method <sup>(2,6,14)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(2,6,13)</sup> 3) Digestion, Flame Atomic Absorption Spectrometric Method <sup>(7,14)</sup> 4) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
7	Chlordane	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(2,9,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup>
8	Chromium	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method <sup>(2,6,14)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(2,6,13)</sup>

3) Digestion,...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
9	Chromium (III)	3) Digestion, Flame Atomic Absorption Spectrometric Method <sup>(7,14)</sup> 4) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup> 1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method; Waste Extraction, Colorimetric Method; Calculation <sup>(2,6,14,16)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method; Waste Extraction, Colorimetric Method; Calculation <sup>(2,6,13,16)</sup> 3) Digestion, Flame Atomic Absorption Spectrometric Method; Alkaline Digestion, Colorimetric Method; Calculation <sup>(7,8,14,16)</sup> 4) Digestion, Inductively Coupled Plasma Method; Alkaline Digestion, Colorimetric Method; Calculation <sup>(7,8,13,16)</sup>
10	Chromium (VI)	1) Waste Extraction, Colorimetric Method <sup>(2,16)</sup> 2) Alkaline Digestion, Colorimetric Method <sup>(8,16)</sup>
11	Cobalt	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(2,6,13)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
12	Copper	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method <sup>(2,6,14)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(2,6,13)</sup> 3) Digestion, Flame Atomic Absorption Spectrometric Method <sup>(7,14)</sup> 4) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
13	2,4-D	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(2,9,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup>
14	DDD	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(2,9,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup>

15 DDE...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
15	DDE	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(2,9,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup>
16	DDT	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(2,9,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup>
17	Dieldrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(2,9,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup>
18	Endrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(2,9,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup>
19	Heptachlor	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(2,9,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup>
20	Lead	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method <sup>(2,6,14)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(2,6,13)</sup> 3) Digestion, Flame Atomic Absorption Spectrometric Method <sup>(7,14)</sup> 4) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
21	Lindane	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(2,9,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup>
22	Mercury	1) Waste Extraction, Digestion, Cold-Vapor Atomic Absorption Spectrometric Method <sup>(2,17)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(2,6,13)</sup>

CONSULTANT COMPANY LIMITED

3) Digestion,...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
		3) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method <sup>(2,18)</sup> 4) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup> 5) Thermal Decomposition Amalgamation and Atomic Absorption Spectrometric Method <sup>(19)</sup>
23	Methoxychlor	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(2,9,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup>
24	Molybdenum	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(2,6,13)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
25	Nickel	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method <sup>(2,6,14)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(2,6,13)</sup> 3) Digestion, Flame Atomic Absorption Spectrometric Method <sup>(7,14)</sup> 4) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
26	Polychlorinated Biphenyls - Aroclor 1016 - Aroclor 1221 - Aroclor 1232 - Aroclor 1242 - Aroclor 1248 - Aroclor 1254 - Aroclor 1260 - 2-Chlorobiphenyl - 2,3-Dichlorobiphenyl - 2,2',5-Trichlorobiphenyl - 2,4',5-Trichlorobiphenyl - 2,2',3,5'-Tetrachlorobiphenyl - 2,2',5,5'-Tetrachlorobiphenyl - 2,3',4,4'-Tetrachlorobiphenyl - 2,2',3,4,5'-Pentachlorobiphenyl	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(2,9,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup>

UN

GC

- 2,2',4,5,5'...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
	- 2,2',4,5,5'-Pentachlorobiphenyl - 2,3,3',4',6-Pentachlorobiphenyl - 2,2',3,4,4',5'-Hexachlorobiphenyl - 2,2',3,4,5,5'-Hexachlorobiphenyl - 2,2',3,5,5',6-Hexachlorobiphenyl - 2,2',4,4',5,5'-Hexachlorobiphenyl - 2,2',3,3',4,4',5-Heptachlorobiphenyl - 2,2',3,4,4',5,5'-Heptachlorobiphenyl - 2,2',3,4,4',5,6-Heptachlorobiphenyl - 2,2',3,4',5,5',6-Heptachlorobiphenyl - 2,2',3,3',4,4',5,5',6-Nonachlorobiphenyl Pentachlorophenol	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(2,9,28)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup> Electrometric Method <sup>(31,32)</sup>
28	pH	
29	Selenium	1) Waste Extraction, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>(2,6,20)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(2,6,13)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup> 4) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>

CONSULTANT COMPANY LIMITED

30 Silver...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
30	Silver	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(2,6,13)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
31	Thallium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(2,6,13)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
32	Toxaphene	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(2,9,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup>
33	Trichloroethylene	1) Waste Extraction, Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(2,12,25)</sup> 2) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
34	Vanadium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(2,6,13)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
35	Zinc	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method <sup>(2,6,14)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(2,6,13)</sup> 3) Digestion, Flame Atomic Absorption Spectrometric Method <sup>(7,14)</sup> 4) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>

ดิน จำนวน 125 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Acenaphthene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,24)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
2	Acetone	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>

UN

GC

3 Aldrin...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
3	Aldrin	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
4	Anthracene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,26)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
5	Antimony	Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
6	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>(7,13)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
7	Atrazine	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
8	Barium	Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
9	Benz(a)anthracene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,26)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
10	Benzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
11	Benzo(b)fluoranthene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,26)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
12	Benzo(k)fluoranthene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,26)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
13	Benzoic acid	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
14	Benzo(a)pyrene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,26)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>

15 Benzo(g,h,i)perylene...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
15	Benzo(g,h,i)perylene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
16	Beryllium	Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
17	Bis(2-chloroethyl)ether	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
18	Bis(2-ethylhexyl)phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
19	Bromodichloromethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
20	Bromoform	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
21	Butanol	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
22	Butyl benzyl phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
23	Cadmium	1) Digestion, Flame Atomic Absorption Spectrometric Method <sup>(7,14)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
24	Carbazole	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
25	Carbon disulfide	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
26	Carbon tetrachloride	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
27	Chlordane	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
28	p-Chloroaniline	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
29	Chlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
30	Chlorodibromomethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>

31 Chloroform...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
31	Chloroform	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
32	2-Chlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
33	Chromium	1) Digestion, Flame Atomic Absorption Spectrometric Method <sup>(7,14)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
34	Chromium (III)	1) Digestion, Flame Atomic Absorption Spectrometric Method; Alkaline Digestion, Colorimetric Method; Calculation <sup>(7,8,14,16)</sup> 2) Digestion, Inductively Coupled Plasma Method; Alkaline Digestion, Colorimetric Method; Calculation <sup>(7,8,13,16)</sup>
35	Chromium (VI)	Alkaline Digestion, Colorimetric Method <sup>(8,16)</sup>
36	Chrysene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,26)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
37	Cyanide	Extraction, Distillation, Colorimetric Method <sup>(28,29,30)</sup>
38	2,4-D	Ultrasonic Extraction, Gas Chromatographic Method <sup>(27)</sup>
39	DDD	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
40	DDE	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
41	DDT	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
42	Dibenz(a,h)anthracene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,26)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>

43 Di-n-butyl phthalate...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
43	Di-n-butyl phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
44	1,2-Dichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
45	1,3-Dichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
46	1,4-Dichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
47	3,3'-Dichlorobenzidine	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
48	1,1-Dichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
49	1,2-Dichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
50	1,1-Dichloroethylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
51	cis-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
52	trans-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
53	2,4-Dichlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
54	1,2-Dichloropopane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
55	1,3-Dichloropropane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
56	1,3-Dichloropropene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
57	Dieldrin	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
58	Diethyl phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
59	2,4-Dimethylphenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>

60 2,4-Dinitrophenol...



ลำดับ	สารมลพิษ	วิธีวิเคราะห์
60	2,4-Dinitrophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
61	2,4-Dinitrotoluene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
62	2,6-Dinitrotoluene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
63	Di-n-Octyl phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
64	Endosulfan	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
65	Endrin	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
66	Ethylbenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
67	Fluoranthene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,24)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
68	Fluorene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,24)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
69	Heptachlor	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
70	Heptachlor epoxide	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>

71 Hexachlorobenzene...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
71	Hexachlorobenzene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
72	Hexachloro-1,3-butadiene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
73	n-Hexane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
74	$\alpha$ -HCH	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
75	$\beta$ -HCH	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
76	$\gamma$ -HCH	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
77	Hexachlorocyclopentadiene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
78	Hexachloroethane	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
79	Indeno(1,2,3-cd)pyrene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,24)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
80	Isophorone	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
81	Lead	1) Digestion, Flame Atomic Absorption Spectrometric Method <sup>(7,14)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
82	Manganese	1) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>

83 Mercury...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
83	Mercury	1) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method <sup>(18)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup> 3) Thermal Decomposition Amalgamation and Atomic Absorption Spectrometric Method <sup>(18)</sup>
84	Methanol	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
85	Methoxychlor	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
86	Methyl bromide	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
87	Methylene chloride	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
88	2-Methylphenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
89	2-Methylnaphthalene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
90	Methyl tert-butyl ether	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
91	Naphthalene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,24)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
92	Nickel	1) Digestion, Flame Atomic Absorption Spectrometric Method <sup>(7,14)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
93	Nitrobenzene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
94	N-Nitrosodiphenylamine	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
95	N-Nitrosodi-n-propylamine	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>

96 Polychlorinated Biphenyls...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
96	Polychlorinated Biphenyls - Aroclor 1016 - Aroclor 1221 - Aroclor 1232 - Aroclor 1242 - Aroclor 1248 - Aroclor 1254 - Aroclor 1260 Polychlorinated Biphenyls - 2-Chlorobiphenyl - 2,3-Dichlorobiphenyl - 2,2',5-Trichlorobiphenyl - 2,4',5-Trichlorobiphenyl - 2,2',3,5'-Tetrachlorobiphenyl - 2,2',5,5'-Tetrachlorobiphenyl - 2,3',4,4'-Tetrachlorobiphenyl - 2,2',3,4,5'- Pentachlorobiphenyl - 2,2',4,5,5'- Pentachlorobiphenyl - 2,3,3',4,6- Pentachlorobiphenyl - 2,2',3,4,4',5'- Hexachlorobiphenyl - 2,2',3,4,5,5'- Hexachlorobiphenyl - 2,2',3,5,5',6- Hexachlorobiphenyl - 2,2',4,4',5,5'- Hexachlorobiphenyl - 2,2',3,3',4,4',5- Heptachlorobiphenyl - 2,2',3,4,4',5,5'- Heptachlorobiphenyl - 2,2',3,4,4',5',6- Heptachlorobiphenyl	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,23)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>  Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,23)</sup>

- 2,2',3,4',5,5',6...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
97	- 2,2',3,4',5,5',6-Heptachlorobiphenyl - 2,2',3,3',4,4',5,5',6-Nonachlorobiphenyl Pentachlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
98	Phenanthrene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,24)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
99	Phenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
100	Pyrene	1) Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,24)</sup> 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
101	Selenium	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method <sup>(7,22)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
102	Silver	Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
103	Styrene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
104	1,1,2,2-Tetrachloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
105	Tetrachloroethylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
106	Toluene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
107	Toxaphene	Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup>
108	TPH (C <sub>9</sub> -C <sub>18</sub> )	1) Purge and Trap, Gas Chromatographic Method <sup>(12,21)</sup> 2) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
109	TPH (C <sub>9</sub> -C <sub>16</sub> )	Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup>
110	TPH (C <sub>16</sub> -C <sub>33</sub> )	Ultrasonic Extraction, Gas Chromatographic Method <sup>(10,22)</sup>
111	1,2,4-Trichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>

112 1,1,1-Trichloroethane...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
112	1,1,1-Trichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
113	1,1,2-Trichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
114	Trichloroethylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
115	2,4,5-Trichlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
116	2,4,6-Trichlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup>
117	1,3,5-Trimethylbenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
118	Vanadium	Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>
119	Vinyl acetate	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
120	Vinyl chloride	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
121	m-Xylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
122	o-Xylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
123	p-Xylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
124	Xylene (Total)	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(12,25)</sup>
125	Zinc	1) Digestion, Flame Atomic Absorption Spectrometric Method <sup>(7,14)</sup> 2) Digestion, Inductively Coupled Plasma Method <sup>(7,13)</sup>

## เอกสารอ้างอิง

- กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2549. เรื่อง กำหนดค่าปริมาณเขม่าควันที่เจือปนในอากาศที่ระบายออกจากปล่องของหม้อไอน้ำที่ใช้เชื้อเพลิงแข็ง. 4 ธันวาคม 2549. เล่มที่ 123 ตอนพิเศษ 113 ก.ราชกิจจานุเบกษา.
- กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2549. เรื่อง กำหนดค่าปริมาณเขม่าควันที่เจือปนในอากาศที่ระบายออกจากปล่องของหม้อไอน้ำที่ใช้เชื้อเพลิงแข็ง. 4 ธันวาคม 2549. เล่มที่ 123 ตอนพิเศษ 113 ก.ราชกิจจานุเบกษา.

3. สมาคมวิศวกรรม...

- สมาคมวิศวกรรมสิ่งแวดล้อมแห่งประเทศไทย. คู่มือวิเคราะห์น้ำเสีย. พิมพ์ครั้งที่ 4. กรุงเทพมหานคร: 2547.
- APHA, AWWA, WEF. Standard Methods for the Examination of Water and Wastewater, 23<sup>rd</sup> ed. Washington, DC: APHA, 2017.
- United States Environmental Protection Agency. Standards of Performance for New Stationary Sources. 40 CFR 60. Appendix A, 2019.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. SW-846, 1997.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Acid Digestion of Sediments, Sludges, and Soils. SW-846 Method 3050B, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Alkaline Digestion for Hexavalent Chromium. SW-846 Method 3060A, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste 3. Physical/Chemical Methods. Separatory Funnel Liquid-Liquid Extraction. SW-846 Method 3510C, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Ultrasonic Extraction. SW-846 Method 3550C, 2007.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Purge and Trap for Aqueous Samples. SW-846 Method 5030C, 2003.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Closed System Purge and Trap and Extraction for Volatile Organics in Soil and Waste Sample. SW-846 Method 5035A, 2000.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Inductively Coupled Plasma-Optical Emission Spectrometry. SW-846 Method 6010D, 2014.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Flame Atomic Absorption Spectrophotometry. SW-846 Method 7000B, 2007.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Arsenic (Atomic Absorption Gaseous Hydride). SW-846 Method 7061A, 1992.

UNITED ANALYST AND ENVIRONMENTAL  
CONSULTANT COMPANY LIMITED

16. United States...

- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Chromium, Hexavalent (Colorimetric). SW-846 Method 7196A, 1992.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Mercury in Liquid Waste (Manual Cold Vapor Technique). SW-846 Method 7470A, 1994.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Mercury in Solid or Semisolid Waste (Manual Cold-Vapor Technique). SW-846 Method 7471B, 1998.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Mercury in Solids and Solutions by Thermal Decomposition, Amalgamation, and Atomic Absorption Spectrophotometry. SW-846 Method 7473, 2007.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Selenium (Atomic Absorption, Borohydride Reduction). SW-846 Method 7742, 1994.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Nonhalogenated Organics Using GC/FID. SW-846 Method 8015D, 2003.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Organochlorine Pesticides by Gas Chromatography. SW-846 Method 8081B, 2007.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Polychlorinated Biphenyls (PCBs) by Gas Chromatography. SW-846 Method 8082A, 2007.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Polynuclear Aromatic Hydrocarbons. SW-846 Method 8100, 1980.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Volatile Organic Compounds by Gas Chromatography/Mass Spectrometry. SW-846 Method 8260D, 2018.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Semivolatile Organic Compounds by Gas Chromatography/Mass Spectrometry. SW-846 Method 8270E, 2018.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Chlorinated Hydrocarbons by Gas Chromatography/Mass Spectrometry. SW-846 Method 8270F, 2018.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Pentafluorobenzoylation Derivatization. SW-846 Method 8270G, 2018.

UNITED ANALYST AND ENVIRONMENTAL  
CONSULTANT COMPANY LIMITED

28. United States...

28. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. **Total and Amenable Cyanide : Distillation**. SW-846 Method 9010C, 2004.

29. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. **Cyanide Extraction Procedure for Solids and Oils**. SW-846 Method 9013A, 2014.

30. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. **Cyanide in Waters and Extracts using Titrimetric and Manual Spectrophotometric Procedures**. SW-846 Method 9014, 2014.

31. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. **pH Electrometric Measurement**. SW-846 Method 9040C, 2004.

32. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. **Soil and Waste pH**. SW-846 Method 9045D, 2004.





