

ภาคผนวก จ

ใบแสดงการตรวจเทียบเครื่องมือ



Request Service No. 101/66

Page 1 of 3

Calibration Certificate

Nomenclature : Brand : Mettler Toledo Type : Top-Loading Electronic Balance

Model : AB204-S Serial No. : 1123163292 (209359)

Submitted by : Laboratory of SECOT CO., LTD.

Location of Calibration : BAL Room , 6th Floor, Secot Co., Ltd.

Calibration range : 0 – 200 g Scale division : 0.0001 g (220 g)

Calibration date : May 25, 2023

Reference Standard No. M220177, M2302167S, M2303005N

Traceable to : Metrological Center SCI ECO Services Company Limited.

: THAI CALIBRATION SERVICES CO., LTD.

Ambient Condition : Temperature 25.20 - 25.50 °C

Humidity 49.7 – 53.4 % RH

Calibrated By : Janista Kui-on Approved By : Narisa Poowasanpetch

(Miss Janista Kui-on)

(Miss Narisa Poowasanpetch)

Testing Officer

Chief of Technical Management

Date : 25/05/2023

Date : 25/05/2023

Issued Date : May 25, 2023

Measurement Report

Request Service No. 101/66

Page 2 of 3

Description : Brand : Mettler Toledo

Type : Top-Loading Electronic Balance

Model : AB204-S

Serial No. : 1123163292 (209359)

Calibration range : 0 – 200 g

Scale division : 0.0001 g (220 g)

Calibration date : May 25, 2023

Ambient Condition : Temperature 25.20-25.50 °C Relative humidity 49.7-53.4 % RH

Measurement data :

1. Repeatability of Reading :

Load (g)	Standard Deviation of Reading (g)	Maximum Difference between Successive Reading (g)
50	0.00006	0.0002
100	0.00008	0.0002
150	0.00010	0.0003
200	0.00007	0.0002

2. Off-Center Loading :

A Mass of 50.0000 g was placed and moved to various position on the pan.

Unit : g

Center	Front	Left	Back	Right	Center	Maximum Difference
50.00002	50.00004	50.00022	50.00012	49.99998	49.99998	0.00024

Issued Date : May 25, 2023

3. Departure from Nominal Value :

Reading (g)	Correction (g)	Uncertainty (+/- g)
0	0.00000	± 0.00008
1	-0.00004	± 0.00008
5	-0.00003	± 0.00012
10	+0.00002	± 0.00010
20	-0.00007	± 0.00009
40	-0.00010	± 0.00012
60	-0.00022	± 0.00016
80	-0.00037	± 0.00014
100	-0.00007	± 0.00016
120	-0.00083	± 0.00018
140	-0.00078	± 0.00021
160	-0.00084	± 0.00023
180	-0.00085	± 0.00026
200	-0.00082	± 0.00027

Calibrated by : Janista Kui-on

(Miss Janista Kui-on)

Testing Officer

Date : 25/05/2023

Approved By : Miss Narisa Poowasanpetch

(Miss Narisa Poowasanpetch)

Chief of Technical Management

Date : 25/05/2023

Issued Date : May 25, 2023



Calibration Certificate

Nomenclature : Brand : Sartorius Type : Top-Loading Electronic Balance

Model : BP 2100 Serial No. : 71003651

Submitted by : Laboratory of SECOT CO., LTD.

Location of Calibration : BAL Room , 6th Floor, Secot Co., Ltd.

Calibration range : 0 - 2000 g Scale division : 0.1 g

Calibration date : January 16, 2023

Reference Standard No. M220177 , M210183 , 22M1064

Traceable to : Metrological Center SCIECO Services Company Limited , Technology Promotion Association (Thailand-Japan)

Ambient Condition : Temperature 25.00-25.50 °C

Humidity 50.9-53.9 % RH

Calibrated By : Khemchuda Insorn

(Miss Khemchuda Insorn)

Testing Officer

Date : 17/01/2023

Approved By : Miss Narisa Poowasanpetch

(Miss Narisa Poowasanpetch)

Chief of Technical Management

Date : 17/01/2023

Issued Date : Jan 17, 2023

Measurement Report

Request Service No. 001/66

Page 2 of 3

Description: Brand : Sartorius

Type : Top-Loading Electronic Balance

Model : BP2100

Serial No. : 71003651

Calibration range : 0 – 2000 g

Scale division : 0.1 g

Calibration date : January 16, 2023

Ambient Condition : Temperature 25.00-25.50°C Relative humidity 50.9-53.9 % RH

Measurement data :

1. Repeatability of Reading :

Load (g)	Standard Deviation of Reading (g)	Maximum Difference between Successive Reading (g)
500	0.042	0.10
1000	0.047	0.20
1500	0.032	0.10
2000	0.032	0.10

2. Off-Center Loading :

A Mass of 1000.00 g was placed and moved to various position on the pan.

Unit : g

Center	Front	Left	Back	Right	Center	Maximum Difference
999.58	999.58	999.62	999.56	999.52	999.56	0.060

Issued Date : Jan 17, 2023

Request Service No.001/66

Page 3 of 3

3. Departure from Nominal Value :

Reading (g)	Correction (g)	Uncertainty (+/- g)
0	0.00	± 0.08
100	0.04	± 0.08
200	0.08	± 0.07
500	0.24	± 0.08
700	0.28	± 0.07
1000	0.32	± 0.07
1200	0.38	± 0.07
1500	0.58	± 0.07
1700	0.56	± 0.08
2000	0.82	± 0.07

Calibrated by :

Khemchuda Insorn

Approved By :

Narisa Poowasanpetch

(Miss Khemchuda Insorn)

(Miss Narisa Poowasanpetch)

Testing Officer

Chief of Technical Management

Date : 17/01/2023

Date : 17/01/2023

Issued Date : Jan 17, 2023



BECTHAI BANGKOK EQUIPMENT & CHEMICAL CO., LTD.
CALIBRATION LABORATORY

300 Phaholyothin Road, Phayathai, Bangkok 10400, Thailand Tel: +66 2615-2929 Fax: +66 2615-2350-1
 E-mail: bkk@becthai.com Website: www.becthai.com



Certificate No. : CAL-23-150

Page : 1 of 4

CERTIFICATE OF CALIBRATION

Equipment : Spectrophotometer
 Manufacturer : Thermo Scientific
 Model : Genesys 150 UV-VIS
 Serial No. : 9A5Y332022
 ID No. : N/A
 Customer : Secot Company Limited
 : 239 Rimklongprapa Road,
 : Bangsue, Bangkok 10800, Thailand
 Location : Laboratory Room
 Date of Receipt : 27 February 2023
 Date of Calibration : 27 February 2023
 Date of Issue : 8 March 2023
 Ambient Temperature : (25±10) °C
 Relative Humidity : (60±20) %
 Condition As-Received : Used Item

Calibrated by

Anusit Boonmee

(Mr. Anusit Boonmee)

Calibration Engineer

Approved by

Jintana Sangthaijaroenlap

(Ms. Jintana Sangthaijaroenlap)

Calibration Manager

The reported expanded uncertainty of measurement was based on a combined standard uncertainty multiplied by a coverage factor $k=2$, providing a level of confidence of approximately 95%.

This certificate may not be reproduced other than in full, except with the prior written approval of the head of Calibration Laboratory.

Indicated values are valid for the state of the Spectrophotometer at the time of calibration only.



BECTHAI BANGKOK EQUIPMENT & CHEMICAL CO., LTD.
CALIBRATION LABORATORY

300 Phaholyothin Road, Phayathai, Bangkok 10400, Thailand Tel: +66 2615-2929 Fax: +66 2615-2350-1
 E-mail: bkk@becthai.com Website: www.becthai.com



Certificate No. : CAL-23-150

Page : 2 of 4

CALIBRATION REPORT

Conditions of this result of calibration

1. Reference Standard Material :

<u>Material</u>	<u>Model</u>	<u>Serial No.</u>	<u>Cert.No.</u>	<u>Due date</u>
Holmium Glass Filter	RM-HG	12705	98236	12 Feb 24
Didymium Glass Filter	RM-DG	13498	98233	12 Feb 24
Neutral Density Filter	RM-1N2N3N	8323	98259	13 Feb 24
Potassium Dichromate Solution	RM-06	23429	98252	12 Feb 24

2. Traceability : This certification is traceable to the International System of Unit maintained at;

The Starna Scientific Ltd. Accredited Calibration Laboratory No. 0659.

3. Method of calibration :

The calibration procedure was carried out according to ASTM E275-08 (2022) and ASTM E925-09 (2014).

4. Result of calibration :

(✓) without adjustment

() after adjustment

5. Equipment Specifications:

Spectral Bandwidth :	2	nm
Data Interval :	0.2	nm
Scan Speed :	Slow	nm/min



BECTHAI BANGKOK EQUIPMENT & CHEMICAL CO., LTD.
CALIBRATION LABORATORY

300 Phaholyothin Road, Phayathai, Bangkok 10400, Thailand Tel: +66 2615-2929 Fax: +66 2615-2350-1
 E-mail: bkk@becthai.com Website: www.becthai.com



Certificate No. : CAL-23-150

Page : 3 of 4

CALIBRATION REPORT

Wavelength Calibration

Certified Values of Reference Material (nm)	Nominal Value (nm)	UUC* Reading (nm)	Error (nm)	Uncertainty of Measurement (\pm nm)
241.74	241.74	241.955	0.215	0.16
637.98	637.98	637.751	-0.229	0.17
879.27	879.27	879.075	-0.195	0.16

Photometric Calibration for Visible

Wavelength (nm)	Certified Values of Reference Material (A)	UUC* Reading (A)	Error (A)	Uncertainty of Measurement (\pm A)
420.0	Zero	0.000	0.0000	0.0028
	0.5716	0.573	0.0014	0.0044
	0.7358	0.733	-0.0028	0.0040
	1.0713	1.073	0.0017	0.0039
440.0	Zero	0.000	0.0000	0.0028
	0.561	0.562	0.0010	0.0042
	0.718	0.715	-0.0030	0.0037
	1.0459	1.047	0.0011	0.0037
465.0	Zero	0.000	0.0000	0.0028
	0.5111	0.512	0.0009	0.0044
	0.6618	0.660	-0.0018	0.0035
	0.9635	0.965	0.0015	0.0034
546.1	Zero	0.000	0.0000	0.0028
	0.5222	0.523	0.0008	0.0036
	0.6687	0.667	-0.0017	0.0031
	0.9768	0.978	0.0012	0.0043
590.0	Zero	0.000	0.0000	0.0028
	0.5541	0.554	-0.0001	0.0035
	0.6975	0.695	-0.0025	0.0031
	1.0206	1.021	0.0004	0.0044
635.0	Zero	0.000	0.0000	0.0028
	0.5398	0.540	0.0002	0.0035
	0.6658	0.664	-0.0018	0.0033
	0.9741	0.974	-0.0001	0.0044

Remark : Each individual filter is measured against the empty filter holder (blank) used to zero the Spectrophotometer.

Note:

UUC* : Unit Under Calibration



BECTHAI BANGKOK EQUIPMENT & CHEMICAL CO., LTD.
CALIBRATION LABORATORY

300 Phaholyothin Road, Phayathai, Bangkok 10400, Thailand Tel: +66 2615-2929 Fax: +66 2615-2350-1
 E-mail: bkk@becthai.com Website: www.becthai.com



Certificate No. : CAL-23-150

Page : 4 of 4

CALIBRATION REPORT

Photometric Calibration for UV

Wavelength (nm)	Certified Values of Reference Material (A)	UUC* Reading (A)	Error (A)	Uncertainty of Measurement (\pm A)
235.0	Zero	0.000	0.0000	0.0050
	0.7345	0.735	0.0005	0.0075
257.0	Zero	0.000	0.0000	0.0050
	0.8498	0.849	-0.0008	0.0074
313.0	Zero	0.000	0.0000	0.0050
	0.2853	0.286	0.0007	0.0055
350.0	Zero	0.000	0.0000	0.0050
	0.6306	0.629	-0.0016	0.0063

Remark : The Potassium Dichromate Filled cells are measured against a Perchloric acid blank.

Note:

UUC* : Unit Under Calibration

- End of Report -

Calibration Certificate

Certificate No.: 2304081-003-01
Client name: SECOT CO., LTD.
Address: 239 Rimklongprapa Road,
 Bangsue, Bangsue, Bangkok 10800

Page 1 of 3

Equipment: CHAMBER (Hot Air Oven)

Manufacturer: BINDER

Model: ED 53

Serial No.: 01-27152

ID No.: N/A


Order No.: 2304081

Operation No.: 2304081-003

Date of Receipt: 27 July 2023

Date of Calibration: 27 July 2023

Calibrated by Mr. Worapob Sooktong
 Scientist

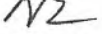Approved by 
 (Mr. Pheraphat Tuanjit)
 Manager, Division of Calibration Laboratory
 Responsible for the Technical Management Team

Date of Issue: 7 August 2023

The uncertainties are for a confidence probability of approximately 95 %.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65



nfi.or.th

Calibration Report

Certificate No.: 2304081-003-01
Equipment: CHAMBER (Hot Air Oven)
 Model: ED 53 Serial No.: 01-27152
 Resolution: 1 °C ID No.: N/A
 Manufacturer: BINDER

Date of Calibration: 27 July 2023

Page 2 of 3

Location: Laboratory, SECOT CO., LTD.

Environment Condition:
 Ambient Temperature (32 ± 1) °C
 Relative Humidity (52 ± 2) %
 Line Voltage (228 ± 1) Volt

Condition of this results of Calibration:

- This instrument was calibrated by insert 9 standard thermometer into its chamber and calibration according to W-TE-014 Based on TLAS G-20-1/02-08 (E): Guidelines for Calibration and Checks of Temperature Controlled Enclosures.
 - The temperature scale used was based on ITS - 90.
 - All data show below were final values and the initial data may be obtained upon request.

2. Reference Standard Instrument :

Instrument	Model	Serial No./ID No.	Certificate No.	Due Date	Through
Digital Thermometer with sensor	34972A	MY49016894	TE 660380-01	22 April 2024	NATIONAL FOOD INSTITUTE
	RTD	CH#101-109/ RTD#101-109			

- This certificate is traceable to International System of Units (SI Units).
- This certificate was certified only for the instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.
- Condition of Calibrated item : Good

UUC Description :

Time of Record 1 Hour 9 Minute At 104, 110 and 180 °C
 Fresh air Damper ☐ Open Position ☐
☒ Close
☐ Not Available

- Result of Calibration : ☒ Without adjustment ☐ After adjustment

F-CS-012 Revision: 01 Date: 20-04-65



nfi.or.th

Calibration Report

Certificate No.: 2304081-003-01
Equipment: CHAMBER (Hot Air Oven)
Model: ED 53 **Serial No.:** 01-27152
Resolution: 1 °C **ID No.:** N/A
Manufacturer: BINDER

Date of Calibration: 27 July 2023
Calibration point: 104, 110 and 180 °C

Page 3 of 3

Calibration result:

Calibration Condition	Temperature (°C)	Relative Humidity (%)	Line Voltage (Volt)
MIN	31.7	50.3	227.1
MAX	32.7	53.5	228.5

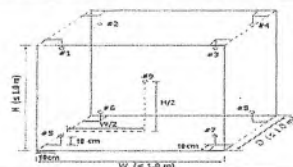


Table1 : Reporting of Temperature

Calibration point (°C)	Measured Temperature (°C) @ Sensor No. (Sensor No.9 is REF)									Uncertainty ± (°C)
	# 1	# 2	# 3	# 4	# 5	# 6	# 7	# 8	# 9	
104	104.79	105.05	104.60	104.30	104.35	103.88	104.29	103.87	103.82	0.78
110	111.06	111.10	110.65	110.38	110.01	109.70	109.80	109.76	109.80	0.80
180	181.06	181.08	180.58	180.53	180.43	180.25	179.97	180.71	180.08	0.90

Table 2 : Reporting of Characterization Result

UUC* Setting (°C)	UUC* reading (°C)			Stability ± (°C)	Uniformity (°C)	Overall Variation (°C)
	MIN	MAX	Average			
104	104	104	104	0.22	1.23	1.55
110	110	110	110	0.25	1.30	1.80
177	177	177	177	0.32	0.99	1.54

Note The quoted uncertainty include " Stability " and " Loading effect (20% of Temp Uniformity) "

UUC* = Unit Under Calibration

Stability = One-half of the greatest maximum difference of measured temperatures at any one sensors, for at least half an hour after reaching steady state.

Uniformity = The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time.

Overall Variation = The difference of the maximum and minimum measured temperatures throughout observation time.

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor $k=2$, providing a level of confidence of approximately 95 %.

----- End -----

Handwritten signature



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
53/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000-27 FAX. 0-2719-9484



Cert.No.: 23CH4
Page.: 1 of 3

Certificate of Calibration

Equipment : pH Meter
Manufacturer : Hanna
Model : HI98190
Serial No. : 06470022101
ID No. : pH No.19
Condition As-Received: Used Item
Received Date : 03 January 2023
Calibration Date : 04 January 2023
Reference : 2301-0006DN-1
Submitted by : Secot Co.,Ltd.
239 Rimklongprapa Road,
Bangsue, Bangkok 10800

Ambient Temperature : (25 ± 2.5) °C
Relative Humidity : (50 ± 15) %
Calibration Procedure : In - house method :
- CP-CH5 by direct measurement with standard voltage calibrator and direct measurement with certified reference material (CRM)
- CP-CH8 by comparison with standard thermometer

Calibrated by : Warakorn Lernagatrakul

Approved by : *Handwritten signature*
Approved Signatory

(✓) Malee Butkruea
() Saitthip Meangmai
() Warakorn Lernagatrakul

Issue Date : 10 January 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Cert.No.: 23CH4
Page.: 2 of 3

Condition of this calibration result

1. Reference Standard Instrument :-

Instrument	Serial No.	ID No.	Cert. No.	Due Date
1) Ref. Standard Thermometer	4982054	110RC044	2211306	27 Oct 2023

This certification is traceable to the International System of Unit maintained at:-

- Traceable to National Institute of Metrology (Thailand), NIMT

2. Certified Reference Materials : The measurement results are traceable to SI through CPA chem Ltd., ANSI-ASQ National Accreditation Board, Accredited No. AR-1835

Buffer Solution	Manufacturer	Lot No.	Exp. date
pH 4.008	CPA chem	826588	09 July 2024
pH 6.987	CPA chem	823322	20 June 2023
pH 10.008	CPA chem	826590	09 July 2023

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration Results

Function : pH Measurement

Performing three buffers standard curve by using buffer nominal pH (4,7,10)

Unit Under Calibration	Standard pH Buffer Solution	Actual pH Reading	Actual mV Reading (mV)	Uncertainty of pH measurement (±)	Coverage factor k
pH Electrode S/N.: 0920044N	4.008	4.010	157.9	0.0044	2.00
	6.987	6.990	-1.6	0.0086	2.00
	10.008	10.007	-163.7	0.0065	2.00

Remark - Can not connect the BNC because the plug does not match with the socket.



Cert.No.: 23CH4
Page.: 3 of 3

Calibration Results

Function : Temperature Measurement

(*) Without adjustment

This equipment was connected with Temperature Probe;

- Model : HI12963

- Serial No. : 0920044N

Dimension of probe;

- Length : 105 mm.

- Diameter : 14 mm.

- Immersion Depth : 100 mm.

Calibration Point (°C)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty of measurement (± °C)	Coverage factor k
20.0	20.002	20.0	-0.002	0.13	2.00
25.0	25.003	25.0	-0.003	0.13	2.00
30.0	30.005	30.0	-0.005	0.13	2.00
35.0	35.002	35.0	-0.002	0.13	2.00

Remark : - UUC* = Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k , providing a level of confidence of approximately 95 %.

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3 : EQUIPMENT CALIBRATION AND TESTING SERVICES

534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250

TEL. 0-2717-3000-29 FAX. 0-2719-9484

Certificate of Calibration

Cert.No.: 23CH656

Page.: 1 of 2

Equipment : Turbidity Meter
Manufacturer : HANNA
Model : HI 98703
Serial No. : 03090023991
ID. No. : -
Condition As-Received: Used Item
Received Date : 19 May 2023
Calibration Date : 23 May 2023
Reference : 2305-0663DN-1
Submitted by : Secot Co.,Ltd.
239 Rimklongprapa Road,
Bangsue, Bangkok 10800
Ambient Temperature : (25 ± 2.5) °C
Relative Humidity : (50 ± 20) %
Calibration Procedure : In - house method : CP-CH11
based on direct measurement by
using Formazin standard solution
Calibrated by : Walalak Sirithean
Approved by : 
Approved Signatory
() Malee Butkruea
(✓) Saithip Meangmai
() Warakorn Lerrigagtrakul
Issue Date : 25 May 2023

The Uncertainties are for a confidence probability of approximately 95%.

This certificate may not be reproduced other than in full, except with the prior written
approval of the head of Calibration and Testing Equipment Services.



Cert.No.: 23CH656

Page.: 2 of 2

Condition of this calibration result

1. Reference Standard Instruments :

This certification is traceable to the International System of unit (SI unit) through
Technology Promotion Association (Thailand-Japan).

Instruments	Serial No.	ID No.	Certificate No.	Due date
1) Thermo-Hygrograph	1103328	130EC010	22H1313	12 June 2023
2) Electronic Balance	N03679	140RC001	22MM49	20 Sep 2023

2. Standard Material : The Formazin suspension has been prepared gravimetric from

Material	Manufacturer	Lot No.	Assay
1) Hexamethylenetetramine	HIMEDIA	0000493947	99.65%
2) Hydrazinium Sulfate	HIMEDIA	0000522014	99.40%

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration result

Performing four - Formazin suspension standard curve by using 0,20,100,800 NTU
Turbidity Meter Serial Number : 03090023991

Standard Formazine suspension (NTU)	UUC* Reading (NTU)	Uncertainty of Measurement (± NTU)	Coverage Factor k
20	19.9	0.39	2.00
40	39.5	0.40	2.00
100	99.4	0.71	2.00
400	396	1.5	2.05

Remark

- UUC* = Unit Under Calibration
- NTU = Nephelometric Turbidity Units

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage
factor k , providing a level of confidence of approximately 95 %.

-o0o-


Calibration Certificate

Certificate No.: 2304081-001-01
Client name: SECOT CO., LTD.
Address: 239 Rimklongprapa Road,
 Bangsue, Bangsue, Bangkok 10800

Page 1 of 3

Equipment: CHAMBER (Incubator)
Manufacturer: MEMMERT
Model: ICP 400
Serial No.: K406.0004
ID No.: N/A
Order No.: 2304081
Operation No.: 2304081-001
Date of Receipt: 27 July 2023
Date of Calibration: 27 July 2023

Calibrated by Mr. Worapob Sooktong
 Scientist

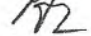Approved by 
 (Mr. Pheraphat Tuanjit)
 Manager, Division of Calibration Laboratory
 Responsible for the Technical Management Team

Date of Issue: 7 August 2023

The uncertainties are for a confidence probability of approximately 95 %.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65



nfi.or.th

Calibration Report

Certificate No.: 2304081-001-01
Equipment: CHAMBER (Incubator)
 Model: ICP 400 Serial No.: K406.0004
 Resolution: 0.1 °C ID No.: N/A
 Manufacturer: MEMMERT
Date of Calibration: 27 July 2023

Page 2 of 3

Location: Laboratory, SECOT CO., LTD.
Environment Condition: Ambient Temperature (30 ± 1) °C
 Relative Humidity (54 ± 1) %
 Line Voltage (228 ± 0) Volt

Condition of this results of Calibration:

- This instrument was calibrated by insert 9 standard thermometer into its chamber and calibration according to W-TE-014 Based on TLAS G-20-1/02-08 (E): Guidelines for Calibration and Checks of Temperature Controlled Enclosures.
 - The temperature scale used was based on ITS - 90.
 - All data show below were final values and the initial data may be obtained upon request.
- Reference Standard Instrument :

Instrument	Model	Serial No./ID No.	Certificate No.	Due Date	Through
Digital Thermometer with sensor	34972A	MY49016894	TE 660380-01	22 April 2024	NATIONAL FOOD INSTITUTE
	RTD	CH#301-309/ RTD#301-309			

- This certificate is traceable to International System of Units (SI Units).
- This certificate was certified only for the instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.
- Condition of Calibrated item : Good

UUC Description :

Time of Record 1 Hour 9 Minute At 20.0 °C
 Fresh air Damper ☐ Open Position ☐
☒ Close
☐ Not Available

- Result of Calibration : ☒ Without adjustment ☐ After adjustment

F-CS-012 Revision: 01 Date: 20-04-65



nfi.or.th

Calibration Report

Certificate No.: 2304081-001-01
Equipment: CHAMBER (Incubator)
Model: ICP 400 **Serial No.:** K406.0004
Resolution: 0.1 °C **ID No.:** N/A
Manufacturer: MEMMERT

Date of Calibration: 27 July 2023

Page 3 of 3

Calibration point: 20.0 °C

Calibration result:

Calibration Condition	Temperature (°C)	Relative Humidity (%)	Line Voltage (Volt)
MIN	28.6	53.0	227.3
MAX	31.4	54.1	228.1

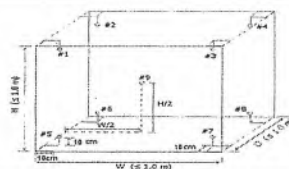


Table1 : Reporting of Temperature

Calibration point (°C)	Measured Temperature (°C) @ Sensor No. (Sensor No.9 is REF)									Uncertainty ± (°C)
	# 1	# 2	# 3	# 4	# 5	# 6	# 7	# 8	# 9	
20.0	20.42	20.39	20.40	20.43	20.47	20.49	20.42	20.41	20.43	0.27

Table 2 : Reporting of Characterization Result

UUC* Setting (°C)	UUC* reading (°C)			Stability ± (°C)	Uniformity (°C)	Overall Variation (°C)
	MIN	MAX	Average			
20.0	20.0	20.1	20.0	0.065	0.053	0.220

Note The quoted uncertainty include " Stability " and " Loading effect (20% of Temp Uniformity) "

UUC* = Unit Under Calibration

Stability = One-half of the greatest maximum difference of measured temperatures at any one sensors, for at least half an hour after reaching steady state.

Uniformity = The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time.

Overall Variation = The difference of the maximum and minimum measured temperatures throughout observation time.

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor $k=2$, providing a level of confidence of approximately 95 %.

----- End -----

Handwritten signature



Calibration Certificate

Certificate No.: 2304081-002-01
Client name: SECOT CO., LTD.
Address: 239 Rimklongprapa Road,
Bangsue, Bangsue, Bangkok 10800

Page 1 of 3

Equipment: Water Bath

Manufacturer: MEMMERT

Model: WB 29

Serial No.: I698.0051

ID No.: N/A

Order No.: 2304081

Operation No.: 2304081-002

Date of Receipt: 27 July 2023

Date of Calibration: 27 July 2023

Calibrated by Mr.Worapob Sooktong
Scientist

Approved by *Handwritten signature*
(Mr.Pheraphat Tuanjit)
Manager, Division of Calibration Laboratory
Responsible for the Technical Management Team

Date of Issue: 7 August 2023

The uncertainties are for a confidence probability of approximately 95 %.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.



Calibration Report

Certificate No.: 2304081-002-01
Equipment: Water Bath
Model: WB 29 Serial No.: I698.0051
Resolution: 0.1 °C ID No.: N/A
Manufacturer: MEMMERT
Date of Calibration: 27 July 2023

Page 2 of 3

Location: Laboratory, SECOT CO., LTD.
Environment Condition: Ambient Temperature (24 ± 1) °C
Relative Humidity (58 ± 2) %
Line Voltage (229 ± 1) Volt

Condition of this results of Calibration:

- This instrument was calibrated by insert 5 standard thermometer into its liquid bath and calibration according to W-TE-011 based on ASTM E715-80 (2016): Standard Specification for Gravity-Convection and Forced-Circulation Water Baths.
- The temperature scale used is ITS - 90.
- All data show below were final values and the initial data may be obtained upon request.

2. Reference Standard Instrument :

Instrument	Model	Serial No./ID No.	Certificate No.	Due Date	Through
Digital Thermometer with sensor	34972A	MY49016894	TE 660380-01	22 April 2024	NATIONAL FOOD INSTITUTE
	RTD	RTD#201-205 / CH#201-205			

- This certificate is traceable to International System of Units (SI Units).
- This certificate was certified only for the instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.
- Condition of Calibrated item : Good

UUC Description:

Time of Record 1 Hour 9 Minute At 95.0 °C

7. Result of Calibration :
- | | |
|---|--------------------|
| X | Without adjustment |
| | After adjustment |



nfi.or.th


Calibration Report

Certificate No.: 2304081-002-01
Equipment: Water Bath
Model: WB 29 Serial No.: I698.0051
Resolution: 0.1 °C ID No.: N/A
Manufacturer: MEMMERT
Date of Calibration: 27 July 2023

Page 3 of 3

Calibration point: 95.0 °C
Calibration result:

Calibration Condition	Temperature (°C)	Relative Humidity (%)	Line Voltage (Volt)
Min	23.0	56.3	227.5
Max	25.0	60.2	229.6

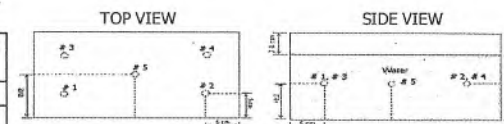

Sensor Installation Location

Table1 : Reporting of Temperature

Calibration Point (°C)	Measured Temperature (°C) @ Sensor No. (Sensor No.5 is REF)					Uncertainty ± (°C)
	# 1	# 2	# 3	# 4	# 5	
95.0	95.03	94.96	95.10	94.97	95.02	0.28

Table 2 : Reporting of Characterization Result

UUC* Setting (°C)	UUC* reading (°C)			Stability ± (°C)	Uniformity (°C)	Overall Variation (°C)
	MIN	MAX	Average			
95.0	94.9	95.1	95.0	0.18	0.080	0.47

Note The quoted uncertainty include " Stability " and " Loading effect (20% of Temp Uniformity)"

UUC* = Unit Under Calibration

Stability = One-half of the greatest maximum difference of measured temperatures at any one sensors, for at least half an hour after reaching steady state.

Uniformity = The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time.

Overall Variation = The difference of the maximum and minimum measured temperatures throughout observation time.

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k= 2, providing a level of confidence of approximately 95 %.

----- End -----



nfi.or.th



MAINTENANCE REPORT
ATOMIC ABSORPTION SPECTROPHOTOMETER MODEL
3110 + HGA600 + FIAS100 + AMALGAM

Customer : บริษัท ซีคอต จำกัด Date Tested: June 15, 2566
 Address : 239 ถนนริมคลองประปา Recommendation Recertification
 แขวงบางซื่อ เขตบางซื่อ Period 6 Months
 กรุงเทพฯ 10800 Recertification Due: December 15, 2566
 User Name: คุณ อารยา Date Last Certified: December 16, 2565
 Phone: 02-9593600 ext. 507 Visit Number: 1 OF 2
 E-mail: labmail@secot.co.th TH Onesource Phone: 081-7316733
 E-mail: thonesource@gmail.com

CONFIGURATION TESTED

MODEL	SERIAL NUMBER	SOFTWARE
AA-3110	311N6062102	AAWINLAB 3.2
HGA 600	2698	
AS 60	2124	
FIAS 100	1114	
AMALGAM	160S2110102	

TEST STANDARD USED	PART NUMBER
Copper	N9300183
GFAAS Mixed STD	N9300244
PE standard of Mercury	N9300174



MAINTENANCE REPORT
ATOMIC ABSORPTION SPECTROPHOTOMETER MODEL
3110 + HGA600 + FIAS100 + AMALGAM

SERIAL NUMBER 311N6062102 DATE TESTED June 15, 2566

1. OPTIC CHECKS

A. Optical alignment condition (if necessary) ☐ OK

B. Condition of Mirrors, Lenses etc. ☐ OK

C. D2 and HCL beam adjust (if necessary) ☐ OK

2. ELECTRONICS CHECKS

A. Power Supplies

+ 5.00 Vdc ± 0.2 Vdc	+ 5.0 Vdc
+ 11.50 Vdc ± 0.2 Vdc	+ 11.4 Vdc
+ 15.00 Vdc ± 1.0 Vdc	+ 15.2 Vdc
- 15.00 Vdc ± 1.0 Vdc	- 14.9 Vdc

B. D2 Power supplies

+150 Vdc	NA Vdc
+ 450 Vdc	NA Vdc

C. PMT Power supply

- 250 Vdc	-249.2 Vdc
-----------	------------

3. GAS SYSTEM CHECKS

A. Leak test all internal and external gas box joints ☐ OK

B. All gas box safety features ☐ OK

C. Burner system including nebulizer and all o-ring and gasket ☐ OK

D. Drain system ☐ OK

4. FIAS CHECK

A. Output power supplies

+5 VDC ± 0.25 VDC.	5.01 VDC.	+40 VDC. ± 0.5 VDC.	40.01 VDC.
--------------------	-----------	---------------------	------------

B. Valve and pump clean ☐ OK



MAINTENANCE REPORT
ATOMIC ABSORPTION SPECTROPHOTOMETER MODEL
3110 + HGA600 + FIAS100 + AMALGAM

SERIAL NUMBER 311N6062102DATE TESTED June 15, 2566**5. PERFORMANCE TEST FOR FLAME**

A. Performance Tests with PE standard.

A1. Run Std. Of Cu and Cr at 324.8 ; 357.9 nm, Concentration 4 , 4 ppm respectivelyResults = 0.228, 0.244 Abs, with flow spoiler. respectivelyCharacteristic Concentration 0.077 ; 0.072 mg/L respectivelyA2. Run Std. of Pb at 283.3 nm; Concentration 20 ppmResults = 0.205 Abs, with flow spoiler.Characteristic Concentration 0.429 mg/LB. Performance Tests (For C₂H₂ + N₂O Flame)Run Std. Of Al at 309.3 nm; Concentration 50 ppmResults = 0.217 Abs, with flow spoiler.Characteristic Concentration 1.014 mg/L**6. PERFORMANCE TEST FOR FIAS****ACTUAL VALUE**

A. Characteristic mass for Mercury

(500 ul of 10 ug/l Hg for 0.07 Abs.)

Characteristic Mass 314 pg / 0.0044 Abs.RSD ≤ 2%0.074 Abs.297.3 pg/0.0044 Abs.1.24 %

B. Characteristic mass for Arsenic

(500 ul of 10 ug/l As for 0.45 Abs.)

Characteristic Mass 48 pg / 0.0044 Abs.RSD ≤ 2%0.482 Abs.45.6 pg/0.0044 Abs.1.54 %

C. Characteristic mass for Mercury Amalgamation

(1000 ul of 1.0 ug/l Hg for 0.03 Abs.)

Characteristic Mass 147 pg / 0.0044 Abs.RSD ≤ 2%0.032 Abs.137.5 pg/0.0044 Abs.1.76 %

MAINTENANCE REPORT
ATOMIC ABSORPTION SPECTROPHOTOMETER MODEL
3110 + HGA600 + FIAS100 + AMALGAM

SERIAL NUMBER 311N6062102DATE TESTED June 15, 2566**7. PERFORMANCE CHECK FOR FURNACE**

A. Internal & External gas flow

☐ OK

B. Contract Cylinder (replace if necessary)

☐ OK

C. Quartz Windows

☐ OK

D. Gas Tubing and Joins

☐ OK

E. Cooling System

☐ OK**8. AUTOSAMPLER CHECK**

A. Arm and gears

☐ OK

B. Sample and Rinse Pump

☐ OK

C. Tray and Sensors

☐ OK**9. PERFORMANCE TEST FOR FURNACE****ACTUAL VALUE****Test run using Chromium**1. Standard Deviation after 5 replicates of blank ≤ 0.005 0.004

2. Characteristic mass (5 ug / L for Cr, 3 pg/0.0044 A-s)

3.0 pg / 0.0044 A-s

Peak Area

0.148

A-s

Relative Standard Deviation ≤ 2.0 %0.24

%

Test run using Lead

Characteristic mass (20 ug / L for Pb, 10 pg/0.0044 A-s)

9.3

pg / 0.0044 A-s

Peak Area

0.188

A-s

Relative Standard Deviation ≤ 2.0 %0.64

%



FSR1206

MAINTENANCE REPORT
ATOMIC ABSORPTION SPECTROPHOTOMETER MODEL
3110 + HGA600 + FIAS100 + AMALGAM

SERIAL NUMBER 311N6062102DATE TESTED June 15, 2566

Remarks :

NA Mean no applicant

This is to certify that the above tests have been performed and the configuration tested



meets



does not meet

the PerkinElmer Specifications listed on this certificate.

This certificate does not modify PerkinElmer's standard terms and condition of sale,
including warranty terms.**TH ONE SOURCE CO., LTD.***Krungchai T.*(**Krungchai Treevichien**)

Customer Support Engineer

Page 5 of 5

**Certificate of Calibration****ICS-1000 : Cation , Anion (ID#057)**

This certificate is to verify that instrument below are calibrated

by Archemica Lab Co.,Ltd.

ICS-1000 S/N : 04090295

RFC-30 S/N : 04100666

For**SECOT Co.,Ltd.**Operator Signature : *Amontep Janparn*

Date : Dec 16, 2022

(Mr.Amontep Janparn)

Applications Chemist

Agilent CrossLab Start Up Services

Agilent 5100 5110 ICP-OES Preventive Maintenance

Agilent Preventive Maintenance provides factory recommended service for your analytical instruments to assure reliable operation and the accuracy of your results

Delivered by highly trained and certified service engineers using genuine Agilent parts and supplies, Agilent Preventive Maintenance provides what you need to reduce unplanned downtime and keep your systems operating at their peak performance.

This checklist is used as a guide for completing the preventive maintenance tasks. A signed copy of this checklist is provided for your records.

Introduction

Customer Information

- Customers should provide all necessary operating supplies upon request of the engineer.
- A customer representative should be available to the engineer while performing the preventive maintenance procedures. Customers are responsible for regular maintenance and are encouraged to observe the service representative.
- Any parts not included in the Parts Lists section of this document are not part of the recommended Preventive Maintenance service nor are they included in the price of this service.
- If a system requires the use of extra or special procedures and/or parts for the maintenance service, then these must be ordered separately and charged as a repair, which may incur additional costs.
- For customers using HF applications, the instrument should be returned to its standard sample introduction system.

Important Customer Web Links

- To access **Agilent University**, visit <http://www.agilent.com/crosslab/university/> to learn about training options, which include online, classroom and onsite delivery. A training specialist can work directly with you to help determine your best options.
- To access the **Agilent Resource Center** web page, visit <https://www.agilent.com/en-us/agilentresources>. The following information topics are available:
 - Sample Prep and Containment
 - Chemical Standards
 - Analysis
 - Service and Support
 - Application Workflows
- The **Agilent Community** is an excellent place to get answers, collaborate with others about applications and Agilent products, and find in-depth documents and videos relevant to Agilent technologies. Visit <https://community.agilent.com/welcome>
- Videos about specific preparation requirements for your instrument can be found by searching the **Agilent YouTube** channel at <https://www.youtube.com/user/agilent>
- **Need to place a service call?** Flexible Repair Options | Agilent

Service Engineer's Responsibilities

- Contact the customer and ensure that all necessary supplies are available before the preventive maintenance visit.
- Only select those pages that relate to the system or module being serviced.
- Complete empty fields with the relevant information.
- Complete the relevant checkboxes in the checklist using either a "X" or tick mark "✓".
- Check "**Service not applicable**" check boxes to indicate services/tasks not delivered, as appropriate.
- Complete the Preventive Maintenance services in the most logical order relevant to the individual system service in the order of the tasks listed.
- Complete the **Service Review** section together with the customer.
- Complete the fields for page numbers at the foot of each selected page
- Add relevant page numbers to selected pages and complete the total number of pages field in the Service Completion section
- **Ask the customer to sign the Service Verification section including the customer's and your signature.**

Instrument Maintenance

System Information

- ☐ Check this box if an instrument configuration report is attached instead of completing the table.

Instrument System Name and ID	5110 VDV ICP-OES
Instrument System Site and Location	Secot Co., Ltd.

List System Component Product Numbers	List the Serial Numbers of each Component
1. G 2015 A	MY 16230003
2. G 2410 A	AU 16151341
3. G 8481-20000	3B 1641345
4.	
5.	
6.	
7.	
8.	
9.	

ICP-OES Configuration Table	Circle the type or write in the type if other
Nebulizer Type	SeaSpray OneNeb Conikal Other
Spray Chamber	Cyclonic Single Pass Cyclonic Double Pass Other
Torch	Radial Dual View Other
Torch Type	One Piece Semi Demountable Fully Demountable Other
Injector Diameter	2.4mm 0.8mm 1.4mm 0.8mm Other
Injector Material	Quartz Ceramic Other

Preparation

- ☒ Discuss any specific issues with the customer before starting.
- ☒ Review the instrument logbook for recorded problems and comments.
- ☒ Save instrument control settings before starting the procedure.
- ☒ Perform a general inspection of the system for cleanliness.
- ☒ Check for proper installation of parts, assemblies, sensors etc.
- ☒ Check system for required installation of components and implementation of Service Notes
- ☒ Check for required firmware/software updates and verify with customers if they would like them installed.
- ☐ For HF application systems, if standard sample introduction system was not installed, ask the customer to install it. N/A
- ☒ Ask the customer to remove any samples from the ICP-OES sample introduction area, auto sampler or around the ICP-OES.

Preventive Maintenance Procedures

Record Pre-PM instrument performance

- ☒ Run Instrument Performance test.
- ☒ Record results in Instrument Performance Test Results Table – Pre-PM.

Clean and inspect ICP-OES system

- ☒ Look for any obvious external damage or problems.
- ☒ Inspect water cooling hoses, gas lines and power cord for excessive wear or damage.
- ☒ Perform a general internal inspection of the system for excessive dust accumulation, clean if necessary.
- ☒ Inspect sample introduction components and record any required maintenance in the Service Engineer Comments and notify the customer as the required actions required.
- ☒ Record the instrument operating conditions in the ICP-OES Status Results Table.
- ☒ Replace the polychromator purge filter.
- ☒ Replace the radial pre-optics window
- ☒ Replace the axial pre-optics window for SVDV and VDV instruments.
- ☒ Check exhaust flow for the correct positive extraction at the exhaust duct to insure they meet minimum specifications.
- ☒ Replace air inlet dust filter.
- ☐ Replace high capacity air inlet dust filter element if installed. **N/A**
- ☒ Remove and clean instrument water inlet filter.

Agilent Water Recirculator

- ☐ Service not applicable
- ☒ Drain cooling fluid and remove any particles from the chiller reservoir
- ☒ Remove, clean and reinstall water inlet metal mesh filter if present.
- ☒ Re fill with Agilent Cool Clear cooling fluid.
- ☒ Clean the cooling system Air filter and the condenser.

SPS 3 Auto Sampler

- ☒ Service not applicable
- ☐ Power cycle the autosampler and verify successful initialization.
- ☐ Inspect X and Z axis belts for wear. Replace is necessary.
- ☐ Clean X and Z axis slide shafts.
- ☐ Using customer's racks and the Agilent software move the sample probe to the 4 outermost corners and rinse port, ensure that the probe is approximately centered in the vial.

SPS 4 Auto sampler

- ☐ Service not applicable
- ☒ Clean the spill tray, rack location mat, end frames and chassis with a damp soft cloth and diluted mild detergent.
- ☒ Clean the auto sampler cover panels, if cover kit is installed, with domestic window cleaner.
- ☒ Check the X-axis and Z-axis drive belts for cracks, splits, damaged teeth, excessive fraying, color changes or degradation from fumes.
- ☒ Check the X-axis, Theta-axis and Z-axis FFC cables for cracks, incorrect positioning, damaged edges or damaged connectors.
- ☒ Pump Tubing Replacement. Replace peristaltic pump tubing. Replace all tubing that goes from the rinse station to the pump and from the pump to the waste/rinse bottles **only check; passed**
- ☒ Test using customer's tray and move the sample probe to the sample vial 1, wash vial and rinse port and ensure that the probe is centered in the vial. If not use calibration wizard and calibrate the position.

AVS 4, 6, 7 Advanced Valve System

- ☒ Service not applicable
- ☐ Replace valve rotor seal
- ☐ Check fittings for signs of leaks
- ☐ Check tubing including autosampler tubing for kinks or excessive wear
- ☐ Check high flow pump for signs of leaks

ICP-OES adjustment

- ☒ Check position of Zn peak, adjust if required.
- ☒ Check Argon Ratio, adjust to specified value if required.
- ☒ Perform Detector Calibration.
- ☒ Perform Instrument Calibration.

Record Post-PM instrument performance

- ☒ Run Instrument Performance test.
- ☒ Record results in Instrument Performance Test Results Table - Post PM.
- ☒ For systems using ICP Expert version 7.3 and above, run the following Instrument tests
 - ☒ Subsystem Communications Test
 - ☒ Air Flow
 - ☒ Water Flow
 - ☒ Gas Flows
 - ☒ RF Generator
 - ☒ Camera Test
 - ☒ Optics Test
 - ☒ Nebulizer Test
- ☒ Record the result in the Instrument Test Results Table

Restore Instrument

- ☐ For HF applications, ask the customer to reinstall their sample introduction system. **W18**
- ☒ Leave system in an idle state: on and purging.
- ☒ Guidance: If the PM service is performed prior to a qualification service, then use the qualification procedure as a guide for final instrument set up and checkout.

Service Review

- ☒ Attach available reports/printouts of all tests to this documentation.
- ☒ Record the Preventive Maintenance service activity in the customer's records/logbook.
- ☒ Record the PM event in the Smart Alerts logbook, if applicable.
- ☒ Update/reset instrument maintenance counters as appropriate.
- ☒ Affix the PM sticker to the system or instrument logbook based on the customer's request.
- ☒ Complete the Service Engineer Comments section if there are additional comments.
- ☒ Review this service, parts replaced, and test results obtained with the customer.
- ☒ If the instrument firmware was updated, record the details of the change in the Service Engineer's Comments box. Systems in a compliant environment may need additional documentation.
- ☒ Complete the Signature Page with both Service Engineer and Customer signatures.

Test Results

Instrument Performance Test Results Table

Note: These measurements do not form part of any specification and are for reference only.

	Pre PM Sensitivity Check		Post PM Sensitivity Check	
	Radial	Axial *	Radial	Axial*
Zn 213.857 nm SRBR	5017.6	12139.3	4952.3	12610.2
Mn 257.610 nm SRBR	12915.3	38934.0	12860.9	41076.8
Al 396.152 nm SBR	12.5	29.0	12.0	25.4
K 766.491 nm SBR	5.1	33.6	4.3	62.6

* Axial result is not applicable for G8016AA, G8012AA Radial View instruments.

Instrument Test Results Table

Note: The Instrument Test results are for systems using ICP Expert version 7.3 and above only.

Instrument Test	Result
Subsystem Communications Test	Pass
Air Flow	Pass
Water Flow	Pass
Gas Flows	Pass
RF Generator	Pass
Camera Test	Pass
Optics Test	Pass
Nebulizer test	Pass

ICP-OES Status Results Table

Note: These measurements do not form part of any specification and are for reference only.

Measurement	Standby Mode		Plasma On	
Mains Voltage	218.805	VAC	217.170	VAC
Mains Current	0.099	A	0.228	A
Instrument Temperature	23.8	°C	24.2	°C
RF Air Flow (sensor speed)	12.0	Hz	18.0	Hz
Plasma Exhaust Temperature	No measurement		30.6	°C
Water Flow Oscillator	No measurement		1.45	L/min
Water Flow Detector	1.12	L/min	1.09	L/min
Water Inlet Temperature	25.9	°C	21.6	°C
Polychromator Temperature	35.0	°C	35.0	°C
CCD Temperature	-20.0	°C	-39.8	°C
Thermal Stabilizer	35.0	°C	35.0	°C
Argon Supply Pressure	598.26	kPa	557.05	kPa
Purge Gas Supply Pressure*1	595.68	kPa	589.30	kPa
Option Gas Supply Pressure*1	-	kPa	-	kPa
Nebulizer Flow	No measurement		0.70	L/min
Nebulizer Back Pressure	No measurement		310.48	kPa
Plasma Gas Flow	No measurement		11.98	L/min
Auxiliary Gas Flow	No measurement		1.00	L/min
RF Power	No measurement		1199.4	W
RF Supply Current	No measurement		8.655	A
RF Supply Voltage	No measurement		184.774	V

*1 If option installed

Consumed PM Parts

Part Description	Part Number	Product or Model# where used	Quantity consumed
Axial Pre-Optic Window	G8010-68014	G8010A, G8011A, G8014A/G8015A	1
Radial Pre-Optic Window	G8010-68015	All	1
Agilent Cool Clear Coolant Fluid	5799-0037	Agilent Water Recirculator	-
Purge Gas Filter	G8010-60136	All	1
Air Inlet filter	G8000-68002	All	1
High Capacity Air Filter	G8010-60189	Optional	-
Rotor seal for 6-7 port valve for AVS6/7	G8494-60002	G8494A/G8495	-
Rotor seal for 4 port valve for AVS4	G8493-60002	G8493A	-
Rinse solution to rinse station 2.5mm id x 1m	G8410-80123	SPS 4	-
Barb connector 2.5mm-1.5mm ID	G8410-80124	SPS 4	-
PVC waste tubing, 8mm od x 5mm id, 2m	G8410-80122	SPS 4	-
Additional Parts may be required from engineer's stock:			
X axis drive belt	5410047500	SPS 3	-
Z axis drive belt	5410047400	SPS 3	-
Peristaltic pump tubing, PVC SolvaFlex, 3 bridged,	3710049000	SPS 4	-

Consumed Parts Reference

(Purchased by customer, not included as part of PM)

☐ Section Not Applicable.

Part Description	Part Number	Product or Model# where used	Quantity consumed
------------------	-------------	------------------------------	-------------------

Signature Page

Service Engineer Comments (optional)

If there are any specific points you wish to note as part of performing the installation or other items of interest for the customer, please write in this box.

Service Verification

Service Request Number:

6006343519

Service Engineer Name:

Kanyakorn S.

Service Engineer Signature:

Kanyakorn S.

Total number of pages in this document:

14


Date Service Completed:

31 Jul 2023

Customer Name:

Customer Signature:

Arach Tisparuch

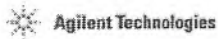

Certificate of Completion

Learner Name:	Kanyakorn Sukpathrajareon
Title Of Course:	ANV-CE-ICPOES-2-008-A: Agilent 5100 ICP-OES Support Neophyte Training
Completion Date:	November 2, 2017
Certified By Company:	Learning at Agilent

All Service and Support training certificates have the following specific limitations.


A certificate for Service and Support training is only valid while employed by Agilent Technologies or while working as an Agilent-authorized service provider, through which the service employee has ongoing access to Agilent's: Safety Alerts, Service Notes, internal technical updates, update training, current documentation, technical support, current parts, and parts updates. Completion of training alone, without being employed by Agilent Technologies, does not qualify an individual to safely install, service or maintain Agilent products.



Agilent CrossLab Start Up Services Agilent 7890 Gas Chromatograph Preventive Maintenance Checklist

Agilent Preventive Maintenance provides factory recommended service for your analytical instruments to assure reliable operation and the accuracy of your results.

Delivered by highly trained and certified service engineers using genuine Agilent parts and supplies, Agilent Preventive Maintenance provides everything you need to reduce unplanned downtime and keep your systems operating at their peak. This checklist will be completed at the end of the service and provided to you as a record of the preventive maintenance activities.

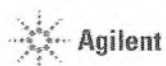

Agilent 7890 GC Preventive Maintenance Checklist



Introduction

Customer Information

- Customers should provide all necessary operating supplies upon request of the engineer.
- A customer representative should be available to the engineer while performing the preventive maintenance procedures.
- Any parts, not included in the Parts Lists section of this document, are not part of the recommended Preventive Maintenance service, nor are they included in the price of this service.
- If a system requires the use of extra or special procedures and/or parts for the maintenance service, then these must be ordered separately and charged as a repair, which may incur additional costs.

Important Customer Web Links

- For more information about **Agilent Technologies services**, please visit our website using the following URL: <http://www.agilent.com/en-us/products/crosslab-instrument-services/service-repair>
- The **Agilent Community** is an excellent place to get answers, collaborate with others about applications and Agilent products, and find in-depth documents and videos relevant to Agilent technologies. Visit <https://community.agilent.com/welcome>.
- To access **Agilent University**, visit <http://www.agilent.com/crosslab/university/> to learn about training options, which include online, classroom and onsite delivery. A training specialist can work directly with you to help determine your best options.
- A useful **Agilent Resource Center** web page is available, which includes short videos on maintenance, quick lists of consumables for new instruments, and other valuable information. Check out the Resource Page here: <https://www.agilent.com/en-us/agilentresources>.
- Need technical support, FAQs, supplies? – visit our **Support Home page** <http://www.agilent.com/search/support>.
- **Videos** about specific preparation requirements for your instrument can be found by searching the **Agilent YouTube** channel at <https://www.youtube.com/user/agilent>.
- **7890B Manuals** are also available on Agilent.com:
 - o **Safety**
https://www.agilent.com/cs/library/usermanuals/public/7890B_Safety.pdf
 - o **Installation and First Startup**
https://www.agilent.com/cs/library/usermanuals/Public/7890B_Installation.pdf
 - o **Operation Manual**
https://www.agilent.com/cs/library/usermanuals/Public/7890B_Operation.pdf
 - o **Maintaining Your GC**
https://www.agilent.com/cs/library/usermanuals/public/G3430-90052%207890B_Maintaining%20Guide.pdf



Service Engineer's Responsibilities

- Contact the customer and ensure that all necessary supplies are available before the preventive maintenance visit.
- Only select those pages that relate to the system or module being serviced.
- Complete empty fields with the relevant information.
- Complete the relevant checkboxes in the checklist using either a "X" or tick mark "✓".
- Check "Section not applicable" check boxes to indicate services/tasks not delivered, as appropriate.
- Complete the Preventive Maintenance service in the order of the tasks listed.
- Complete the Service Review section together with the customer.
- Complete the fields for page numbers at the foot of each selected page
- Complete the total number of pages field in the Service Completion section
- Ask the customer to sign the Service Completion section including the customer's and your signature.

Additional Instruction Notes

- Check for any active service notes for this unit. If there are any applicable "Safety" or "Modification Recommended" Service notes, plan to implement the changes on this unit before doing any qualification service.
- Do not implement firmware updates, unless you get approval from the customer and are sure that they are compatible with the instrument control software.

System Information

- ☒ Check this box if an instrument configuration report is attached instead of completing the table below.

Instrument System Name and ID

GC7890B

CN15343147

Instrument System Site and Location

Secot Co., Ltd.

Instrument room.

List System Component Product Numbers

List the Serial Numbers of each Component

 1. G3440B
 2. G4513A
 3. G4514A

 CN15343147
 CN1910080
 CN19080006

Preparation

- ☒ Discuss any specific issues with the customer before starting.
- ☒ Review the instrument logbook for recorded problems and comments.
- ☒ Save instrument control settings before starting the procedure.
- ☒ Perform a general inspection of the system for cleanliness.
- ☒ Check for proper installation of parts, assemblies, sensors etc.
- ☒ Check system for required installation of components, settings as defined by current Service Notes.
- ☒ Check for required firmware updates and verify with customers if they would like them installed.
- ☒ Before starting the following procedures, record the Detector Signal Output(s) in the results table. If the GC is turned OFF or in a service mode, comparing the detector outputs before and after the service is not possible.

Preventive Maintenance Procedure

Clean and inspect GC

- ☒ Unplug power cord from the power source.
- ☒ Open GC covers and vacuum/remove any dust/debris. Pay particular attention to cooling fans.
- ☒ Inspect internal connectors for proper contact and placement.
- ☒ Reconnect Power to the GC. Power the GC on and verify the power on self-test passed.
- ☒ Verify oven motor spins freely and turns on with the oven door closed; off when the door is opened.
- ☒ Verify operation of all other fans - the inlet and EPC cooling fans.
- ☒ Verify oven intake/outlet flap assembly is operating smoothly while heating and cooling the oven

Inlet and detector consumable replacement

- ☒ For the inlets installed, perform inlet maintenance as defined in the 7890 manual – "Maintaining Your GC" - for the Inlet(s) installed.
- ☒ Replace the split vent trap cartridge filter on units with these inlets: Split/Splitless Capillary (SSL), Multi-Mode Inlet (MMI), Programmed Temperature Vaporizer (PTV), Volatiles Interface (VI).
- ☒ If the inlet system is used in Split Mode with viscous samples, inspect and clean the split vent tube on the inlet and flush or replace the tubing between the inlet and the split vent trap.
- ☒ If the GC includes a Flame Ionization Detector (FID), replace the jet. If the ignitor shows any buildup of sample or corrosion, replace the ignitor. Examine the FID collector and castle assemblies for contamination – clean as necessary.

Zero Sensors and Leak test

- ☒ Zero all pressure sensors per the procedure in the 7890 "Advanced User Guide".
- ☒ Perform inlet pressure decay test(s) as defined in the 7890 "Troubleshooting Manual".
If the PM is done in preparation for an Operational Qualification, then the pressure decay test defined within that protocol can be used for the PM.
- ☒ Record if test passed or failed in the results table.

ALS Maintenance

- ☐ Section NOT applicable
- ☒ Check all cabling and configuration settings between GC, tray, and injectors.
- ☒ Vacuum or remove any dust, especially around fans.
- ☒ Check operation of all fans.
- ☒ Check syringe for smooth plunger operation.
- ☒ Check for smooth operation of the needle support – clean if necessary

Restore Instrument

- ☒ Restore the normal operating conditions or customer method using the Data System.
- ☒ Purge the system with carrier flow for 15 minutes
- ☒ Bake out the system, then restore the normal operating conditions
- ☒ After equilibration, check and record the post PM detector signal output values.
Results should be similar or lower than the detector outputs recorded prior to PM.
- ☐ Perform a chemical checkout. If this is a routine PM, inject the customer's sample using the ALS if applicable. This will act as a final checkout of both the ALS and the GC.

Note: If the PM Service is performed prior to a qualification service, then use the qualification procedure as a guide for final instrument set up and checkout.

Signature Page

Service Review

- ☐ Attach available reports/printouts of all tests to this documentation.
☒ Record the Preventive Maintenance service activity in the customer's records/logbook.
☒ Update/reset instrument maintenance counters as appropriate.
☒ Affix the PM sticker to the system or instrument logbook based on the customer's request.
☒ Complete the Service Engineer Comments section if there are additional comments.
☒ Review with the customer this service, parts replaced, and test results obtained.
☒ If the instrument firmware was updated, record the details of the change in the Service Engineer's Comments box or if necessary, in the customer's IQ records.
☐ Supply the customer with a copy of the Smart Alerts flyer.
☐ Describe Smart Alerts to the customer.
☐ Install Smart Alerts if requested.

7890 GC Test Results Table

Detector Signal Outputs	Before PM Service	After PM Service
Front detector output <i>AFCD</i>	<i>N/A</i>	<i>126.2</i>
Back detector output <i>PID</i>	<i>N/A</i>	<i>22.6</i>
AUX detector output	<i>N/A</i>	<i>N/A</i>
Pressure decay test	Expected test result	Actual test result
Front inlet pressure decay test	Pass	<i>Pass</i>
Back inlet pressure decay test	Pass	<i>Pass</i>

7890 Parts List Table

The following kits are recommended for capillary and purged packed inlets. If this is a general PM and the customer has a preferred set of consumables, you may use the customer's consumables.

Part description	Part number	Product or model# where used	Quantity consumed
SSL Capillary Inlet PM kit, Splitless	5188-6497	7890A/B	<i>N/A</i>
SSL Capillary Inlet PM kit, split	5188-6496	7890A/B	<i>2</i>
SSL Capillary Ultra Inert Inlet Gold Seal with Washer	5190-6144	7890A/B	<i>N/A</i>
SSL Capillary Ultra Inert Inlet Splitless Liner - Single taper with Glass Wool	5190-2293	7890A/B	<i>N/A</i>
SSL Capillary Ultra Inert Inlet Low Pressure Drop Split Liner - with Glass Wool	5190-2295	7890A/B	<i>N/A</i>
PP Inlet PM kit	5188-6498	7890A/B	<i>N/A</i>
Split vent trap PM kit, single cartridge (for MMI, PTV & VI)	5188-6495	7890A/B	<i>N/A</i>
MMI Cleaning Kit	G3510-60820	7890A/B	<i>N/A</i>
PTV Septumless Head Rebuild Kit	5182-9747	7890A/B	<i>N/A</i>
PTV Septumless Head Teflon Guide	5182-9748	7890A/B	<i>N/A</i>
Ignitor (glow plug) assembly with O-ring	19231-60680	7890A/B	<i>1</i>
FID Collector Rebuild/Cleaning Kit	G1531-67000	7890A/B	<i>N/A</i>
Standard .011-inch FID Jet for capillary FID base	G1531-80560	7890A/B	<i>1</i>
High Temperature .018-inch FID Jet for capillary FID base	G1531-80620	7890A/B	<i>N/A</i>
Standard .018-inch FID Jet for packed column with packed FID base	18710-20119	7890A/B	<i>N/A</i>
Standard .011-inch FID Jet for capillary column with packed/adaptable FID base	19244-80560	7890A/B	<i>N/A</i>
High Temperature .018-inch FID Jet for capillary column with packed/adaptable FID base	19244-80620	7890A/B	<i>N/A</i>
NPD Jet, universal fit, .011-inch ID	G1534-80580	7890A/B	<i>N/A</i>
NPD Jet, universal fit, .011-inch ID Extended tip	G1534-80590	7890A/B	<i>N/A</i>
SSL Capillary Ultra Inert Inlet Gold Seal with Washer	5190-6144	7890A/B	<i>N/A</i>
SSL Capillary Ultra Inert Inlet Splitless Liner - Single taper with Glass Wool	5190-2293	7890A/B	<i>N/A</i>
**FID Collector Replacement Kit, if needed	G1531-67001	7890A/B	<i>N/A</i>

Service Engineer Comments

If there are any specific points you wish to note as part of performing the service or other items of interest for the customer, please write include them in this box.

Service Completion

Service request number

600641153

Agilent signature

SAH

Date service completed


29 May 2023

Customer signature

Mum. Wansapthi

Total number of pages in this document

9 pages

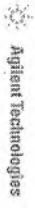

Certificate of Completion

Learner Name:

Saenguthai Saeng Tarek

Title Of Course:

AN-ASP/CE/CSE-GC-1-001-M: 7890/7820 GC and OL GC Standalone Chemstation I&F/Service

Completion Date:

November 23, 2014

Certified By Company:

Learning at Agilent

All Service and Support training certificates have the following specific limitations.

A certificate for Service and Support training is only valid while employed by Agilent Technologies or while working as an Agilent-authorized service provider, through which the service employee has ongoing access to Agilent's Safety Alerts, Service Notes, internal technical updates, update training, current documentation, technical support, current parts, and parts updates. Completion of training alone, without being employed by Agilent Technologies, does not qualify an individual to safely install, service or maintain Agilent products.



Certificate of Completion

Learner Name:	Saenguthai Saeng Tarak
Title Of Course:	AN-CE-GC-IL-022-A: Advanced GC Detectors Application and Troubleshooting Labs
Completion Date:	November 25, 2014
Certified By Company:	Learning at Agilent

All Service and Support training certificates have the following specific limitations.

A certificate for Service and Support training is only valid while employed by Agilent Technologies or while working as an Agilent-authorized service provider, through which the service employee has ongoing access to Agilent's Safety Alerts, Service Notes, internal technical updates, update training, current documentation, technical support, current parts, and parts updates. Completion of training alone, without being employed by Agilent Technologies, does not qualify an individual to safely install, service or maintain Agilent products.

Certificate of Completion

Learner Name:	Saenguthai Tarak
Title Of Course:	AN-CE-GC-IL-022-A: Advanced GC Detectors Application and Troubleshooting
Completion Date:	March 18, 2016
Certified By Company:	Learning at Agilent

All Service and Support training certificates have the following specific limitations.

A certificate for Service and Support training is only valid while employed by Agilent Technologies or while working as an Agilent-authorized service provider, through which the service employee has ongoing access to Agilent's Safety Alerts, Service Notes, internal technical updates, update training, current documentation, technical support, current parts, and parts updates. Completion of training alone, without being employed by Agilent Technologies, does not qualify an individual to safely install, service or maintain Agilent products.


Agilent CrossLab Start Up Services

Agilent 7890 Gas Chromatograph

Preventive Maintenance Checklist

Agilent Preventive Maintenance provides factory recommended service for your analytical instruments to assure reliable operation and the accuracy of your results.

Delivered by highly trained and certified service engineers using genuine Agilent parts and supplies, Agilent Preventive Maintenance provides everything you need to reduce unplanned downtime and keep your systems operating at their peak. This checklist will be completed at the end of the service and provided to you as a record of the preventive maintenance activities.

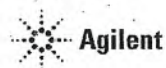

Agilent 7890 GC Preventive Maintenance Checklist


Introduction

Customer Information

- Customers should provide all necessary operating supplies upon request of the engineer.
- A customer representative should be available to the engineer while performing the preventive maintenance procedures.
- Any parts, not included in the Parts Lists section of this document, are not part of the recommended Preventive Maintenance service, nor are they included in the price of this service.
- If a system requires the use of extra or special procedures and/or parts for the maintenance service, then these must be ordered separately and charged as a repair, which may incur additional costs.

Important Customer Web Links

- For more information about **Agilent Technologies services**, please visit our website using the following URL: <http://www.agilent.com/en-us/products/crosslab-instrument-services/service-repair>
- The **Agilent Community** is an excellent place to get answers, collaborate with others about applications and Agilent products, and find in-depth documents and videos relevant to Agilent technologies. Visit <https://community.agilent.com/welcome>.
- To access **Agilent University**, visit <http://www.agilent.com/crosslab/university/> to learn about training options, which include online, classroom and onsite delivery. A training specialist can work directly with you to help determine your best options.
- A useful **Agilent Resource Center** web page is available, which includes short videos on maintenance, quick lists of consumables for new instruments, and other valuable information. Check out the Resource Page here: <https://www.agilent.com/en-us/agilentresources>.
- Need technical support, FAQs, supplies? – visit our **Support Home page** <http://www.agilent.com/search/support>.
- **Videos** about specific preparation requirements for your instrument can be found by searching the **Agilent YouTube** channel at <https://www.youtube.com/user/agilent>.
- **7890B Manuals** are also available on Agilent.com:
 - o **Safety**
https://www.agilent.com/cs/library/usermanuals/public/7890B_Safety.pdf
 - o **Installation and First Startup**
https://www.agilent.com/cs/library/usermanuals/Public/7890B_Installation.pdf
 - o **Operation Manual**
https://www.agilent.com/cs/library/usermanuals/Public/7890B_Operation.pdf
 - o **Maintaining Your GC**
https://www.agilent.com/cs/library/usermanuals/public/G3430-90052%207890B_Maintaining%20Guide.pdf

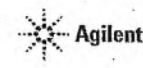

Service Engineer's Responsibilities

- Contact the customer and ensure that all necessary supplies are available before the preventive maintenance visit.
- Only select those pages that relate to the system or module being serviced.
- Complete empty fields with the relevant information.
- Complete the relevant checkboxes in the checklist using either a "X" or tick mark "✓".
- Check "Section not applicable" check boxes to indicate services/tasks not delivered, as appropriate.
- Complete the Preventive Maintenance service in the order of the tasks listed.
- Complete the Service Review section together with the customer.
- Complete the fields for page numbers at the foot of each selected page
- Complete the total number of pages field in the Service Completion section
- **Ask the customer to sign the Service Completion section including the customer's and your signature.**

Additional Instruction Notes

- Check for any active service notes for this unit. If there are any applicable "Safety" or "Modification Recommended" Service notes, plan to implement the changes on this unit before doing any qualification service.
- Do not implement firmware updates, unless you get approval from the customer and are sure that they are compatible with the instrument control software.

System Information

- ☐ Check this box if an instrument configuration report is attached instead of completing the table below.

Instrument System Name and ID	7890A GC System / GCMS
Instrument System Site and Location	SECOT Co., Ltd.

List System Component Product Numbers	List the Serial Numbers of each Component
1. G3440A	CN10750035
2. N/A	N/A
3. N/A	N/A
4. N/A	N/A
5. N/A	N/A
6. N/A	N/A
7. N/A	N/A
8. N/A	N/A
9. N/A	N/A
10. N/A	N/A

Preparation

- ☒ Discuss any specific issues with the customer before starting.
- ☒ Review the instrument logbook for recorded problems and comments.
- ☒ Save instrument control settings before starting the procedure.
- ☒ Perform a general inspection of the system for cleanliness.
- ☒ Check for proper installation of parts, assemblies, sensors etc.
- ☒ Check system for required installation of components, settings as defined by current Service Notes.
- ☒ Check for required firmware updates and verify with customers if they would like them installed.
- ☐ Before starting the following procedures, record the Detector Signal Output(s) in the results table. If the GC is turned OFF or in a service mode, comparing the detector outputs before and after the service is not possible.

Preventive Maintenance Procedure

Clean and inspect GC

- ☒ Unplug power cord from the power source.
- ☒ Open GC covers and vacuum/remove any dust/debris. Pay particular attention to cooling fans.
- ☒ Inspect internal connectors for proper contact and placement.
- ☒ Reconnect Power to the GC. Power the GC on and verify the power on self-test passed.
- ☒ Verify oven motor spins freely and turns on with the oven door closed; off when the door is opened.
- ☒ Verify operation of all other fans - the inlet and EPC cooling fans.
- ☒ Verify oven intake/outlet flap assembly is operating smoothly while heating and cooling the oven

Inlet and detector consumable replacement

- ☒ For the inlets installed, perform inlet maintenance as defined in the 7890 manual – "Maintaining Your GC" - for the inlet(s) installed.
- ☒ Replace the split vent trap cartridge filter on units with these inlets: Split/Splitless Capillary (SSL), Multi-Mode Inlet (MMI), Programmed Temperature Vaporizer (PTV), Volatiles Interface (VI).
- ☒ If the inlet system is used in Split Mode with viscous samples, inspect and clean the split vent tube on the inlet and flush or replace the tubing between the inlet and the split vent trap.
- ☐ If the GC includes a Flame Ionization Detector (FID), replace the jet. If the ignitor shows any buildup of sample or corrosion, replace the ignitor. Examine the FID collector and castle assemblies for contamination – clean as necessary.

Zero Sensors and Leak test

- ☒ Zero all pressure sensors per the procedure in the 7890 "Advanced User Guide".
- ☒ Perform inlet pressure decay test(s) as defined in the 7890 "Troubleshooting Manual".
If the PM is done in preparation for an Operational Qualification, then the pressure decay test defined within that protocol can be used for the PM.
- ☒ Record if test passed or failed in the results table.

ALS Maintenance

- ☒ Section NOT applicable
- ☐ Check all cabling and configuration settings between GC, tray, and injectors.
- ☐ Vacuum or remove any dust, especially around fans.
- ☐ Check operation of all fans.
- ☐ Check syringe for smooth plunger operation.
- ☐ Check for smooth operation of the needle support – clean if necessary

Restore Instrument

- ☒ Restore the normal operating conditions or customer method using the Data System.
- ☒ Purge the system with carrier flow for 15 minutes
- ☒ Bake out the system, then restore the normal operating conditions
- ☒ After equilibration, check and record the post PM detector signal output values.
Results should be similar or lower than the detector outputs recorded prior to PM.
- ☐ Perform a chemical checkout. If this is a routine PM, inject the customer's sample using the ALS if applicable. This will act as a final checkout of both the ALS and the GC.

Note: If the PM Service is performed prior to a qualification service, then use the qualification procedure as a guide for final instrument set up and checkout.

Signature Page

Service Review

- ☒ Attach available reports/printouts of all tests to this documentation.
- ☒ Record the Preventive Maintenance service activity in the customer's records/logbook.
- ☒ Update/reset instrument maintenance counters as appropriate.
- ☒ Affix the PM sticker to the system or instrument logbook based on the customer's request.
- ☒ Complete the Service Engineer Comments section if there are additional comments.
- ☒ Review with the customer this service, parts replaced, and test results obtained.
- ☐ If the instrument firmware was updated, record the details of the change in the Service Engineer's Comments box or if necessary, in the customer's IQ records.
- ☐ Supply the customer with a copy of the Smart Alerts flyer.
- ☐ Describe Smart Alerts to the customer.
- ☐ Install Smart Alerts if requested.

7890 GC Test Results Table

Detector Signal Outputs	Before PM Service	After PM Service
Front detector output	N/A	N/A
Back detector output	N/A	N/A
AUX detector output	N/A	N/A
Pressure decay test	Expected test result	Actual test result
Front inlet pressure decay test	Pass	Pass
Back inlet pressure decay test	Pass	Pass

7890 Parts List Table

The following kits are recommended for capillary and purged packed inlets. If this is a general PM and the customer has a preferred set of consumables, you may use the customer's consumables.



Part description	Part number	Product or model# where used	Quantity consumed
SSL Capillary Inlet PM kit, Splitless	5188-6497	7890A/B	1
SSL Capillary Inlet PM kit, split	5188-6496	7890A/B	1
SSL Capillary Ultra Inert Inlet Gold Seal with Washer	5190-6144	7890A/B	N/A
SSL Capillary Ultra Inert Inlet Splitless Liner - Single taper with Glass Wool	5190-2293	7890A/B	N/A
SSL Capillary Ultra Inert Inlet Low Pressure Drop Split Liner - with Glass Wool	5190-2295	7890A/B	N/A
PP Inlet PM kit	5188-6498	7890A/B	N/A
Split vent trap PM kit, single cartridge (for MMI, PTV & VI)	5188-6495	7890A/B	N/A
MMI Cleaning Kit	G3510-60820	7890A/B	N/A
PTV Septumless Head Rebuild Kit	5182-9747	7890A/B	N/A
PTV Septumless Head Teflon Guide	5182-9748	7890A/B	N/A
Ignitor (glow plug) assembly with O-ring	19231-60680	7890A/B	N/A
FID Collector Rebuild/Cleaning Kit	G1531-67000	7890A/B	N/A
Standard .011-inch FID Jet for capillary FID base	G1531-80560	7890A/B	N/A
High Temperature .018-inch FID Jet for capillary FID base	G1531-80620	7890A/B	N/A
Standard .018-inch FID Jet for packed column with packed FID base	18710-20119	7890A/B	N/A
Standard .011-inch FID Jet for capillary column with packed/adaptable FID base	19244-80560	7890A/B	N/A
High Temperature .018-inch FID Jet for capillary column with packed/adaptable FID base	19244-80620	7890A/B	N/A
NPD Jet, universal fit, .011-inch ID	G1534-80580	7890A/B	N/A
NPD Jet, universal fit, .011-inch ID Extended tip	G1534-80590	7890A/B	N/A
SSL Capillary Ultra Inert Inlet Gold Seal with Washer	5190-6144	7890A/B	N/A
SSL Capillary Ultra Inert Inlet Splitless Liner - Single taper with Glass Wool	5190-2293	7890A/B	N/A
**FID Collector Replacement Kit, if needed	G1531-67001	7890A/B	N/A

Service Engineer Comments

If there are any specific points you wish to note as part of performing the service or other items of interest for the customer, please write include them in this box.

N/A

Service Completion

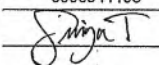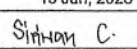Service request number 6006041193 Date service completed 15 Jun, 2023
 Agilent signature  Customer signature 
 Total number of pages in this document 10

Do not include this section/page in the published, customer-facing PDF version.

This page is only relevant for Agilent source documents for document control purposes and is NOT intended for customer viewing. Refer to the SPIIFPM checklist Authoring Guide for more information.

Document Control Logs

Revision Log

Revision	Date	Author	Reason for update
Revision of document	Date of Issuance	Author of document	Author to describe main features/changes made for this specific revision
1.0 Draft	4-Mar-2011	Dave Park	Migrated the content of revision A.01.05 to the new Agilent template. Reviewed by subject matter expert, Dave Park.
1.1 Draft	20-Jan-2015	Dave Park	Added Split Vent trap to MMI, PTV and VE - also PTV and FID PM Parts
1.2 Draft	31-March-2015	Dave Park	Added Ultra Inert Gold Seal and Liner to SS + Consumables
A.01.11	10-Dec-2015	Dave Park	Added step to perform maintenance on the Split Vent Tube and .018" FID Jet part numbers - Fixed broken web links
2.00	30-Dec-2020	Gary Boardman	Updated New Template and terminology change: Familiarization to Introduction. Create New Agile Document Number: D0007063

Approval Log

Revision	Approver	Title of approver
Add revision number	Add approver name here	Add approver's function or title here
A.01.06	Don Gage	Product support manager
A.01.09	Kai Meng	Product support manager
A.01.10	Suneetha Tippireddy	Product support manager
A.01.11	Suneetha Tippireddy	Product support manager
2.00	Josh Roark	GC Product Support Manager

Designated Evaluation Log

Revision	Designated Evaluator (DE)	Title of DE	DE Number
Add revision number	Add name	Add function or title	Add DE number here
2.00	Michael Zumwalt	CrossLab Start Up Services Application Consulting Lead	44166.7597222222



ELECTRICAL AND ELECTRONICS INSTITUTE
FOUNDATION FOR INDUSTRIAL DEVELOPMENT

975 Moo 4, Bangpoo Industrial Estate, Soi 8, Sukhumvit Road km 37,

Phraek Sa, Mueang Samut Prakan, Samut Prakan 10280

Tel: +66 2709 4860 Fax: +66 2324 0917



Certificate No.: CP20230033EA
Operation No.: CP2023010024

Certificate of Calibration

Equipment: Sound Calibrator
Manufacturer: RION
Model/Type: NC-74
Serial No.: 34283648
ID No.: -
Customer: SECOT Co.,Ltd.
Address: 239 Rimklongprapa Rd., Bangsue,
Bangkok 10800 Thailand
Received Date: 10 January 2023
Calibrated Date: 13 January 2023
Issued Date: 16 January 2023
Calibrated by: Ms. Juntaporn Kunhakom

Approved by: _____

(Mr. Sittichai Swaksuriyawong)
Group Manager

This report was prepared electronically using applicable electronic signature. Printing or copy of file are considered as a copy of the document.

The reported uncertainty of measurement was based on standard uncertainty multiplied by a coverage factor (k) providing a level of confidence of approximately 95%. This certificate may not be reproduced other than in full except with the prior written approval of the Electrical and Electronics Institute, Foundation for Industrial Development.



ELECTRICAL AND ELECTRONICS INSTITUTE
FOUNDATION FOR INDUSTRIAL DEVELOPMENT

Certificate No.: CP20230033EA

Calibration Report

Equipment: Sound Calibrator
Manufacturer: RION
Model/Type: NC-74
Serial No.: 34283648
ID No.: -
Ambient Temperature: (23 ± 2) °C
Relative Humidity: (50 ± 15) %
Pressure: (101.3 ± 1.5) kPa

Method of Calibration :-
IEC 60942:2017

Condition of this result of calibration

1. Reference standards instrument :-

Instrument	Model	Serial No.	Cert. No.	Due Date
1) Standard microphone	4180	2661000	AA-1020-22	14 June 2023
2) Waveform Generator	33511B	MY52302264	CK20220058EA	19 June 2023
3) Audio Analyzing DMM	2015-P	4079144	E1U221042	16 March 2023
4) Pressure humidity and Temperature Transmitter	PTU301	F0640002	CL1-P220024 CD20220165EA	17 March 2023 24 July 2023

2. This result of calibration was found accurate as shown on date and place of calibration only.

3. This certification is traceable to the international system of unit maintained at :-

- Reference standards instrument for Acoustic function
- National Institute of Metrology (Thailand)
- Reference standards instrument for Electrical function
- Electrical and Electronics Institute; NSG Accredited Calibration No.0119

Result of Calibration:-

1. Function : Sound pressure level

Norminal Frequency (Hz)	Specified Sound Pressure level (dB)	Measured value (dB)	Deviated value ^[1] (dB)	Acceptance limit ^[3] (dB)
1000	94	94.24	0.24	±0.25

2. Function : Frequency

Norminal Sound Pressure level (dB)	Specified Frequency (Hz)	Measured value (Hz)	Deviated value ^[2] (%)	Acceptance limit ^[3] (%)
94	1000	1003.0	0.3	±0.7



ELECTRICAL AND ELECTRONICS INSTITUTE
FOUNDATION FOR INDUSTRIAL DEVELOPMENT

Certificate No.: CP20230033EA

Calibration Report

3. Function : Total distortion + noise

Nominal Sound Pressure level (dB)	Nominal Frequency (Hz)	Measured value ^[4] (%)	Acceptance limit ^[5] (%)
94	1000	1.3	2.5

Uncertainty of measurement

Function	Uncertainty	Maximum-permitted uncertainty of measurement
Sound pressure level	0.10 dB	0.15 dB
Frequency	0.10 %	0.20 %
Total distortion + noise	0.40 %	0.50 %

- Note: [1] The deviated value is the absolute value of the difference between the measured value and the corresponding specified sound pressure level.
[2] The deviated value is the absolute value of the difference in percent between the measured value and the corresponding specified frequency.
[3] The acceptance limit is for the deviated value.
[4] The measured value is the total distortion + noise, measured over the frequency range from 20 Hz to 20 kHz.
[5] The acceptance limit is for the Measured value.
- Remarks: 1. Using the 1/2-inch microphone adaptor NC-74-002.
2. Acceptance limit was IEC 60942:2017 Class 1.
3. The coverage factor $k = 2.00$

-- End of Report --



ELECTRICAL AND ELECTRONICS INSTITUTE
FOUNDATION FOR INDUSTRIAL DEVELOPMENT

975 Moo 4, Bangpoo Industrial Estate, Soi 8, Sukhumvit Road km 37,

Phraek Sa, Mueang Samut Prakan, Samut Prakan 10280

Tel: +66 2709 4860 Fax: +66 2324 0917



Certificate No.: CP20220368EA

Operation No.: CP2022120011

Certificate of Calibration

Equipment: Sound Calibrator
Manufacturer: Cirrus Research Plc
Model/Type: CR:515
Serial No.: 94296
ID No.: -
Customer: SECOT Co.,Ltd.
Address: 239 Rimklongprapa Rd., Bangsue,
Bangkok 10800 Thailand
Received Date: 14 December 2022
Calibrated Date: 20 December 2022
Issued Date: 23 December 2022
Calibrated by: Ms. Juntaporn Kunhakom

Approved by: _____

(Mr. Sittichai Swaksuriyawong)
Group Manager

This report was prepared electronically using applicable electronic signature. Printing or copy of file are considered as a copy of the document.

The reported uncertainty of measurement was based on standard uncertainty multiplied by a coverage factor (k) providing a level of confidence of approximately 95%. This certificate may not be reproduced other than in full except with the prior written approval of the Electrical and Electronics Institute, Foundation for Industrial Development.

Certificate No.: CP20220368EA

Calibration Report

Equipment: Sound Calibrator
Manufacturer: Cirrus Research Plc
Model/Type: CR:515
Serial No.: 94296
ID No.: -
Ambient Temperature: (23 ± 2) °C
Relative Humidity: (50 ± 15) %
Pressure: (101.3 ± 1.5) kPa

Method of Calibration :-

IEC 60942:2017

Condition of this result of calibration

1. Reference standards instrument :-

Instrument	Model	Serial No.	Cert. No.	Due Date
1) Standard microphone	4180	2661000	AA-1020-22	14 June 2023
2) Waveform Generator	33511B	MY52302264	CK20220058EA	19 June 2023
3) Audio Analyzing DMM	2015-P	4079144	E1U221042	16 March 2023
4) Pressure humidity and Temperature Transmitter	PTU301	F0640002	CL1-P220024 CD20220165EA	17 March 2023 24 July 2023

2. This result of calibration was found accurate as shown on date and place of calibration only.

3. This certification is traceable to the international system of unit maintained at :-

Reference standards instrument for Acoustic function

- National Institute of Metrology (Thailand)

Reference standards instrument for Electrical function

- Electrical and Electronics Institute; NSC Accredited Calibration No.0119

Result of Calibration:-

1. Function : Sound pressure level

Normal	Specified Sound	Measured value	Deviated value ^[1]	Acceptance limit ^[3]
Frequency (Hz)	Pressure level (dB)	(dB)	(dB)	(dB)
1000	94	93.90	-0.10	±0.25

2. Function : Frequency

Normal Sound	Specified Frequency	Measured value	Deviated value ^[2]	Acceptance limit ^[3]
Pressure level (dB)	(Hz)	(Hz)	(%)	(%)
94	1000	1000.3	0.0	±0.7

Certificate No.: CP20220368EA

Calibration Report

3. Function : Total distortion + noise

Normal	Normal	Measured value ^[4]	Acceptance limit ^[5]
Sound Pressure level (dB)	Frequency (Hz)	(%)	(%)
94	1000	0.9	2.5

Uncertainty of measurement

Function	Uncertainty	Maximum-permitted uncertainty of measurement
Sound pressure level	0.10 dB	0.15 dB
Frequency	0.10 %	0.20 %
Total distortion + noise	0.40 %	0.50 %

- Note:
- [1] The deviated value is the absolute value of the difference between the measured value and the corresponding specified sound pressure level.
 - [2] The deviated value is the absolute value of the difference in percent between the measured value and the corresponding specified frequency.
 - [3] The acceptance limit is for the deviated value.
 - [4] The measured value is the total distortion + noise, measured over the frequency range from 20 Hz to 20 kHz.
 - [5] The acceptance limit is for the Measured value.

Remarks: 1. Acceptance limit was IEC 60942:2017 Class 1.

2. The coverage factor $k = 2.00$

-- End of Report --

CERTIFICATE OF CALIBRATION

ISSUED BY Noisemeters


DATE OF ISSUE 16 March 2023 CERTIFICATE NUMBER 189327

NoiseMeters

NoiseMeters
Acoustic House
Bridlington Road
Hunmanby
YO14 0PH
United Kingdom
www.noisemeters.com

Page 1 of 1

Test engineer:
Nigel Smith
Electronically signed:



doseBadge Reader

Instrument

Manufacturer: Cirrus Research plc Serial Number: 95168
Model Number: RC:110A Notes:

Calibration Procedure

The tests were carried out in accordance with the requirements of IEC 60942:2003 where applicable.

Date of Calibration: 16 March 2023

Functionality Results

Function	Result
Keypad	Pass
Battery Power	Pass
Display	Pass
Communication	Pass
2 way IR link	Pass
Clock	Pass

Calibration Results

	Level (dB)	Frequency (Hz)	Distortion (% THD + Noise)
Initial	113.90	999.3	0.61
Adjusted	114.00	999.2	0.61
Uncertainty	± 0.11	± 0.14	± 0.10
Tolerances	± 0.60	± 2.00	± 4.00

Environmental Conditions

Pressure: 99.27 kPa
Temperature: 23.3 °C
Humidity: 37.6 %

Notes

This certificate provides traceability of measurement to the SI system of units and/or to units of measurement realised at the National Physical Laboratory or other recognised national metrology institutes. This certificate may not be reproduced other than in full, except with the prior written approval of the issuing laboratory. The results within this certificate relate only to the items calibrated. The reported expanded uncertainty is based on a standard uncertainty multiplied by a coverage factor $k=2$, providing a coverage probability of approximately 95%.

ภาคผนวก จ

ใบอนุญาตขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
จากกรมโรงงานอุตสาหกรรม



ที่ อก ๐๓๑๐(๑)/ ๑๑ ๐๑ ๖

กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๒๐ กรกฎาคม ๒๕๖๖

เรื่อง ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

เรียน กรรมการผู้จัดการ บริษัท ซีคอต จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๗ เมษายน ๒๕๖๖

- สิ่งที่ส่งมาด้วย ๑. รายชื่อผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๑ แผ่น
๒. รายชื่อเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๑ แผ่น
๓. ขอบข่ายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๓๙ แผ่น

ตามหนังสือที่อ้างถึง บริษัท ซีคอต จำกัด ขอต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ว-๒๓๙ สถานที่ ตั้งเลขที่ ๒๓๙ ถนนริมคลองประปา แขวงบางซื่อ เขตบางซื่อ กรุงเทพมหานคร ต่อกรมโรงงานอุตสาหกรรม นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้บริษัท ซีคอต จำกัด ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน โดยมีองค์ประกอบดังนี้

- ก. ผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๑๐ ราย ตามสิ่งที่ส่งมาด้วย ๑
ข. เจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๓๘ ราย ตามสิ่งที่ส่งมาด้วย ๒
ค. ขอบข่ายสารมลพิษที่ได้รับขึ้นทะเบียนให้วิเคราะห์ในน้ำเสีย น้ำใต้ดิน อากาศเสีย สิ่งปฏิกูล หรือวัสดุที่ไม่ใช้แล้ว และดิน ตามสิ่งที่ส่งมาด้วย ๓

หนังสือฉบับนี้จะหมดอายุในวันที่ ๒ พฤษภาคม ๒๕๖๙ หากประสงค์จะต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ให้ยื่นคำขอต่ออายุพร้อมเอกสารประกอบคำขอต่อกรมโรงงานอุตสาหกรรมภายใน ๓๐ วัน ก่อนวันสิ้นอายุของหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ซึ่งคำขอต่ออายุดังกล่าวขอรับได้ที่กรมโรงงานอุตสาหกรรม ทั้งนี้ สามารถยื่นคำขอผ่านระบบอิเล็กทรอนิกส์ได้ที่หน้าเว็บไซต์กรมโรงงานอุตสาหกรรม

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

(นายประสม ดำรงพงษ์)

กองวิจัยและเตือนภัยมลพิษโรงงาน

ผู้อำนวยการกองวิจัยและเตือนภัยมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม

โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕

โทรสาร ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๕๙

ไปรษณีย์อิเล็กทรอนิกส์ saraban@diw.mail.go.th



"อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว"



สิ่งที่ส่งมาด้วย ๑

เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท ซีคอต จำกัด

เลขทะเบียน ว-๒๓๙

ที่ อก ๐๓๑๐(๑)/ ๑๑ ๐๑ ๖

ลงวันที่ ๒๐ กรกฎาคม ๒๕๖๖

ก. ผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๑๐ ราย

๑) นายขรรชัย เกรียงไกรอุดม

ทะเบียนเลขที่ ว-๒๓๙-ค-๐๐๐๒

๒) นางสมฤดี เกรียงไกรอุดม

ทะเบียนเลขที่ ว-๒๓๙-ค-๐๐๐๓

๓) นางสาวอารยา ทิพริกซ์

ทะเบียนเลขที่ ว-๒๓๙-ค-๐๐๐๔

๔) นางสาวเชมชุตตา อินทร์ศรี

ทะเบียนเลขที่ ว-๒๓๙-ค-๐๐๐๕

๕) นางสาวปริตตา สมใจ

ทะเบียนเลขที่ ว-๒๓๙-ค-๐๐๐๖

๖) นางสาวอริยญา มาตา

ทะเบียนเลขที่ ว-๒๓๙-ค-๐๐๐๗

๗) นางสาวลดาวัลย์ วงศ์เจริญ

ทะเบียนเลขที่ ว-๒๓๙-ค-๐๐๐๘

๘) นางสาวณัฏฐา ฤกษ์วันดี

ทะเบียนเลขที่ ว-๒๓๙-ค-๐๐๐๙

๙) นางสาวนริสา ภูวสรเพ็ชญ์

ทะเบียนเลขที่ ว-๒๓๙-ค-๐๐๑๐

๑๐) นางสาวศิริวรรณ อิมสง่า

ทะเบียนเลขที่ ว-๒๓๙-ค-๐๐๑๑

เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
บริษัท ซีคอฟ จำกัด เลขทะเบียน ๖-๒๓๙๙
ที่ อก ๐๓๑๐(๑)/ ๑๑ ๐ ๑ ๖ ลงวันที่ ๒๐ กรกฎาคม ๒๕๖๖
ข. เจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๓๘ ราย

สิ่งที่ส่งมาด้วย ๒

๑) นางสาวสุดาพร สุนทร	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๐๑
๒) นางสาวสุธาทิพย์ เทียนเตี้ย	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๐๓
๓) นางสาวสุนันทา ศิริพัฒนานนท์	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๐๔
๔) นายบวร ดีชัยยะ	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๐๕
๕) นางสาวเกศรินทร์ วรเดชวิทยา	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๐๖
๖) นายอนันต์ พิฆานัน	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๐๗
๗) นายชิตพล สมประสงค์	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๐๘
๘) นางสาวศศิธร พรหมประเสริฐ	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๐๙
๙) นายศิวะนนท์ กุลวงษ์	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๑๐
๑๐) นางสาวอลิษา คณิธรานนท์	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๑๑
๑๑) นางสาวสิริวรรณ แก้วชิงดวง	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๑๒
๑๒) นางสาวปัทมวรรณ สุวรรณวิโรจน์	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๑๓
๑๓) นางสาวนันทิยา เจริญเชื้อ	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๑๔
๑๔) นายวัชรกานต์ ประมาคะเต	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๑๕
๑๕) นายทอง เสงฆ์กุล	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๑๖
๑๖) นางสาวกฤษณา จันทุม	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๑๗
๑๗) นางสาวพรนภา บุตรธรรม	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๑๘
๑๘) นางสาวอารีณี อางปลิว	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๑๙
๑๙) นายธนโชติ ช่างล้อ	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๒๐
๒๐) นางสาวพัชรา สมานจันทร์	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๒๑
๒๑) นางสาวจุฑาทิพย์ แจ่มเรือน	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๒๒
๒๒) นางสาวณิสดา กุ้ยอ่อน	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๒๓
๒๓) นายกิตติพงศ์ ตะเกียงสุข	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๒๔
๒๔) นายจิรวัฒน์ โคตรคำหาญ	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๒๕
๒๕) นายชนะพล อัครผล	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๒๖
๒๖) นางสาวทิพย์สุดา วรรณการ	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๒๗
๒๗) นายสิทธิชัย สว่างวงศ์ไชย	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๒๘
๒๘) นายพิษณุ สีนามเพ็ง	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๒๙
๒๙) นายรัตนชัย ขอบทำกิจ	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๓๐
๓๐) นายธนวุฒิ ด่วนแสง	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๓๑
๓๑) นายณัฐชัย ไชยโคตร	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๓๒
๓๒) นายณัฐดนัย กฤษณะโสม	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๓๓
๓๓) นายศุภชัย สุขใหม่	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๓๔
๓๔) นายรอมฎอน เหมหมมาต	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๓๕
๓๕) นางสาวสุภาวดี บัวแก้ว	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๓๖
๓๖) นางสาวมาธิยาณี อาแว	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๓๗
๓๗) นางสาววิระยา ปัจฉิมบุรณ์	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๓๘
๓๘) นางสาวศลิษา อินทรีย์	ทะเบียนเลขที่ ๖-๒๓๙๙-จ-๐๐๓๙

ริกา

เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
บริษัท ซีคอฟ จำกัด เลขทะเบียน ๖-๒๓๙๙
ที่ อก ๐๓๑๐(๑)/ ๑๑ ๐ ๑ ๖ ลงวันที่ ๒๐ กรกฎาคม ๒๕๖๖

สิ่งที่ส่งมาด้วย ๓

ขอขยายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๓๕๕ รายการ

น้ำเสีย จำนวน 45 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Aldrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
2	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[4] 2) Digestion, Inductively Coupled Plasma Method ^[4]
3	Barium	1) Digestion, Direct Nitrous Oxide-Acetylene Flame Method ^[4] 2) Digestion, Inductively Coupled Plasma Method ^[4]
4	α-BHC	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
5	β-BHC	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
6	δ-BHC	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
7	γ-BHC	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
8	Biochemical Oxygen Demand	1) 5-Day BOD Test, Azide Modification Method ^[4] 2) 5-Day BOD Test, Membrane Electrode Method ^[4]
9	Cadmium	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Method ^[4]
10	Chemical Oxygen Demand	1) Open Reflux, Titrimetric method ^[4] 2) Closed Reflux, Colorimetric method ^[4] 3) Closed Reflux, Titrimetric Method ^[4]
11	Chlordane	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
12	Chromium	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Method ^[4]
13	Color	ADMI Weighted-Ordinate Spectrophotometric Method ^[4]
14	Copper	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Method ^[4]
15	Cyanide	Distillation, Colorimetric method ^[4]
16	4,4'-DDD	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
17	4,4'-DDE	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
18	4,4'-DDT	2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4] 1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
19	Dieldrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
20	Endosulfan I	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
21	Endosulfan II	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
22	Endosulfan Sulfate	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
23	Endrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
24	Endrin Aldehyde	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
25	Formaldehyde	Distillation, Colorimetric Method ^[3]
26	Free Chlorine	1) Iodometric Method ^[4] 2) DPD Colorimetric Method ^[4]
27	Heptachlor	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
28	Heptachlor epoxide	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
29	Hexavalent Chromium	1) Colorimetric Method ^[4] 2) Extraction, Air-Acetylene Flame Method ^[4]
30	Lead	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Method ^[4]
31	Manganese	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Method ^[4]
32	Mercury	Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[4]
33	Methoxychlor	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
34	Nickel	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] <i>3) Digestion...</i>

3) Digestion...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
		3) Digestion, Inductively Coupled Plasma Method ^[4]
35	Oil & Grease	1) Liquid-Liquid, Partition-Gravimetric Method ^[4] 2) Soxhlet Extraction Method ^[4]
36	pH	Electrometric Method ^[4]
37	Phenols	1) Distillation, Chloroform Extraction Method ^[4] 2) Distillation, Direct Photometric Method ^[4]
38	Selenium	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[4] 2) Digestion, Inductively Coupled Plasma Method ^[4]
39	Sulfide	1) Iodometric method ^[4] 2) Methylene blue method ^[4]
40	Temperature	Laboratory and Field Methods ^[4]
41	Total Dissolved Solids	Dried at 180 °C ^[4]
42	Total Kjeldahl Nitrogen	1) Macro Kjeldahl Method ^[4] 2) Semi-Micro Kjeldahl Method ^[4]
43	Total Suspended Solids	Dried at 103-105 °C ^[4]
44	Trivalent Chromium	1) Digestion, Direct Air-Acetylene Flame Method; Colorimetric Method; Calculation ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method; Colorimetric Method; Calculation ^[4] 3) Digestion, Inductively Coupled Plasma Method; Colorimetric Method; Calculation ^[4]
45	Zinc	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Method ^[4] <i>3) Digestion...</i>

น้ำใต้ดิน จำนวน 125 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Acenaphthene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
2	Acetone	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
3	Aldrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
4	Anthracene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
5	Antimony	Digestion, Inductively Coupled Plasma Spectrometric Method ^[4]
6	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[4] 2) Digestion, Inductively Coupled Plasma Method ^[4]
7	Atrazine	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
8	Barium	1) Digestion, Direct Nitrous Oxide-Acetylene Flame Method ^[4] 2) Digestion, Inductively Coupled Plasma Spectrometric Method ^[4]
9	Benz(a)anthracene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
10	Benzene	Purge and Trap Gas Chromatographic/Mass spectrometric Method ^[4]
11	Benzo(b)fluoranthene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
12	Benzo(k)fluoranthene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4] 3 มด)

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
13	Benzoic acid	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
14	Benzo(a)pyrene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
15	Benzo(g,h,i)perylene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
16	Beryllium	Digestion, Inductively Coupled Plasma Spectrometric Method ^[4]
17	Bis(2-chloroethyl)ether	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
18	Bis(2-ethylhexyl)phthalate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
19	Bromodichloromethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
20	Bromoform	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
21	Butanol	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
22	Butyl benzyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
23	Cadmium	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Spectrometric Method ^[4]
24	Carbazole	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
25	Carbon disulfide	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
26	Carbon tetrachloride	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4] 3 มด)

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
27	Chlordane	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
28	p-Chloroaniline	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
29	Chlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
30	Chlorodibromomethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
31	Chloroform	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
32	2-Chlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
33	Chromium	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Spectrometric Method ^[4]
34	Chromium (III)	1) Digestion, Direct Air-Acetylene Flame Method; Colorimetric Method; Calculation ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method; Colorimetric Method; Calculation ^[4] 3) Digestion, Inductively Coupled Plasma Spectrometric Method; Colorimetric Method; Calculation ^[4]
35	Chromium (VI)	1) Colorimetric Method ^[4] 2) Extraction, Air-Acetylene Flame Method ^[4]
36	Chrysene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4] <i>สมย</i>

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
37	Cyanide	1) Distillation, Titrimetric Method ^[4] 2) Distillation, Colorimetric Method ^[4]
38	2,4-D	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
39	DDD	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
40	DDE	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
41	DDT	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
42	Dibenz(a,h)anthracene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
43	Di-n-butyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
44	1,2-Dichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
45	1,3-Dichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
46	1,4-Dichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
47	3,3'-Dichlorobenzidine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
48	1,1-Dichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
49	1,2-Dichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4] <i>สมย</i>

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
50	1,1-Dichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
51	cis-1,2-Dichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
52	trans-1,2-Dichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
53	2,4-Dichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
54	1,2-Dichloropropane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
55	1,3-Dichloropropane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
56	1,3-Dichloropropene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
57	Dieldrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
58	Diethyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
59	2,4-Dimethylphenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
60	2,4-Dinitrophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
61	2,4-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
62	2,6-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
63	Di-n-Octyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
64	Endosulfan	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] <i>รวม</i>

2) Liquid-Liquid...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
65	Endrin	2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4] 1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
66	Ethylbenzene	2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4] Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
67	Fluoranthene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
68	Fluorene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
69	Heptachlor	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
70	Heptachlor epoxide	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
71	Hexachlorobenzene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
72	Hexachloro-1,3-butadiene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
73	n-Hexane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
74	α-HCH	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
75	β-HCH	1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] <i>รวม</i>


2) Liquid-Liquid...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
76	γ-HCH	2) Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4] 1) Liquid-Liquid Extraction, Gas Chromatographic Method ^[4] 2) Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
77	Hexachlorocyclopentadiene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
78	Hexachloroethane	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
79	Indeno(1,2,3-cd)pyrene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
80	Isophorone	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
81	Lead	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Spectrometric Method ^[4]
82	Manganese	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Spectrometric Method ^[4]
83	Mercury	Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[4]
84	Methanol	Purge and Trap Gas Chromatographic/ Mass spectrometric Method ^[4]
85	Methoxychlor	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
86	Methyl bromide	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]


87 Methylene chloride...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
87	Methylene chloride	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
88	2-Methylphenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
89	2-Methylnaphthalene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
90	Methyl tert-butyl ether	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
91	Naphthalene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
92	Nickel	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Spectrometric Method ^[4]
93	Nitrobenzene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
94	N-Nitrosodiphenylamine	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
95	N-Nitrosodi-n-propylamine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
96	Polychlorinated Biphenyls - PCB-1016 - PCB-1221 - PCB-1232 - PCB-1242 - PCB-1248 - PCB-1254 - PCB-1260	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
97	Pentachlorophenol	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
98	pH	Electrometric method ^[4]

99 Phenanthrene...

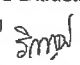ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
99	Phenanthrene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
100	Phenol	1) Distillation, Chloroform Extraction Method ^[4] 2) Distillation, Direct Photometric Method ^[4] 3) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
101	Pyrene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
102	Selenium	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[4] 2) Digestion, Inductively Coupled Plasma Method ^[4]
103	Silver	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Inductively Coupled Plasma Method ^[4]
104	Styrene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
105	1,1,2,2-Tetrachloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
106	Tetrachloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
107	Toluene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
108	TPH (C ₅ -C ₈)	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[12,25]
109	TPH (C ₈ -C ₁₆)	1) Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[9,21] 2) Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass spectrometric Method ^[9,25]
110	TPH (C ₁₆ -C ₃₅)	1) Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[9,21] 

2) Separatory...

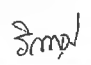ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
		2) Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass spectrometric Method ^[9,25]
111	1,2,4-Trichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
112	1,1,1-Trichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
113	1,1,2-Trichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
114	Trichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
115	2,4,5-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
116	2,4,6-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
117	1,3,5-Trimethylbenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
118	Vanadium	Digestion, Inductively Coupled Plasma Spectrometric Method ^[4]
119	Vinyl acetate	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
120	Vinyl chloride	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
121	m-Xylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
122	o-Xylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
123	p-Xylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4]
124	Xylene (Total)	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ^[4] 

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
125	Zinc	1) Digestion, Direct Air-Acetylene Flame Method ^[4] 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ^[4] 3) Digestion, Inductively Coupled Plasma Spectrometric Method ^[4]

อากาศเสีย (ปล่อยระบาย) จำนวน 27 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Antimony	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
2	Arsenic	1) Isokinetic Sampling, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
3	Beryllium	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
4	Calcium	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
5	Carbon monoxide	Instrumental Analyzer Method ^[5]
6	Chlorine	1) Absorption Sampling, Ion Chromatographic Method ^[5] 2) Isokinetic Sampling, Ion Chromatographic Method ^[5]
7	Chromium	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] <i>3mg/L</i>

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
8	Cobalt	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
9	Copper	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
10	Cresol	Adsorption Sampling, Gas Chromatographic Method ^[5]
11	Dioxin/Furans	Isokinetic Sampling ^[5]
12	Hydrogen chloride	1) Absorption Sampling, Ion Chromatographic Method ^[5] 2) Isokinetic Sampling, Ion Chromatographic Method ^[5]
13	Hydrogen Fluoride	1) Absorption Sampling, Ion Chromatographic Method ^[5] 2) Isokinetic Sampling, Ion Chromatographic Method ^[5]
14	Hydrogen Sulfide	Absorption Sampling, Iodometric Method ^[5]
15	Lead	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
16	Manganese	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
17	Mercury	Isokinetic Sampling, Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[5]
18	Nickel	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] <i>3mg/L</i>

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
19	Opacity	Ringelmann's Method ^[2]
20	Oxides of Nitrogen	1) Absorption Sampling, Phenoldisulfonic acid Method ^[5] 2) Absorption Sampling, Ion Chromatographic Method ^[5] 3) Instrumental Analyzer Method ^[5]
21	Selenium	1) Isokinetic Sampling, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
22	Sulfur dioxide	1) Isokinetic Sampling, Barium-Thorin Titrimetric Method ^[5] 2) Absorption Sampling, Barium-Thorin Titrimetric Method ^[5] 3) Instrumental Analyzer Method ^[5]
23	Sulfuric acid	Isokinetic Sampling, Barium-Thorin Titrimetric Method ^[5]
24	Tin	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
25	Total Suspended Particulate	1) Isokinetic Sampling, Gravimetric Method ^[5] 2) Paired Train, Isokinetic Sampling, Gravimetric Method ^[5]
26	Vanadium	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
27	Xylene	1) Adsorption Sampling, Gas Chromatographic Method ^[5] 2) Adsorption Sampling, Gas Chromatographic/Mass Spectrometric Method ^[5]

สิ่งปฏิกูล...

สิ่งปฏิกูลหรือวัสดุที่ไม่ใช้แล้ว จำนวน 34 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Aldrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,6,9,22] 2) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,6,9,27] 3) Soxhlet Extraction, Gas Chromatographic Method ^[10,22] 4) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27]
2	Antimony	1) Waste Extraction, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 3) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[7,16] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]
3	Arsenic	1) Waste Extraction, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 3) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[7,16] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]
4	Barium	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^[1,6,15]

2) Waste Extraction...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
5	Beryllium	2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 3) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]
6	Cadmium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
7	Chlordane	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^[1,6,15] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 3) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]
8	Chromium	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,9,22] 2) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,27] 3) Soxhlet Extraction, Gas Chromatographic Method ^[10,22] 4) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27] 1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^[1,6,15] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] <i>Simple</i>

3) Digestion...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
9	Chromium (III)	3) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 4) Digestion, Inductively Coupled Plasma Method ^[7,14] 1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method; Waste Extraction, Colorimetric Method; Calculation ^[1,6,15,17] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method; Waste Extraction, Colorimetric Method; Calculation ^[1,6,14,17] 3) Digestion, Flame Atomic Absorption Spectrometric Method; Alkaline Digestion, Colorimetric Method; Calculation ^[7,8,15,17] 4) Digestion, Inductively Coupled Plasma Method; Alkaline Digestion, Colorimetric Method; Calculation ^[7,8,14,17]
10	Chromium (VI)	1) Waste Extraction, Colorimetric Method ^[1,17] 2) Alkaline Digestion, Colorimetric Method ^[8,17]
11	Cobalt	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
12	Copper	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^[1,6,15] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 3) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 4) Digestion, Inductively Coupled Plasma Method ^[7,14] <i>Simple</i>

13 2,4-D...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
13	2,4-D	1) Waste Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,25] 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[25]
14	DDD	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,9,22] 2) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,27] 3) Soxhlet Extraction, Gas Chromatographic Method ^[10,22] 4) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27]
15	DDE	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,9,22] 2) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,27] 3) Soxhlet Extraction, Gas Chromatographic Method ^[10,22] 4) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27]
16	DDT	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,9,22] 2) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,27] 3) Soxhlet Extraction, Gas Chromatographic Method ^[10,22] 4) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27]

17 Dieldrin...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
17	Dieldrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,9,22] 2) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,27] 3) Soxhlet Extraction, Gas Chromatographic Method ^[10,22] 4) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27]
18	Endrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,9,22] 2) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,27] 3) Soxhlet Extraction, Gas Chromatographic Method ^[10,22] 4) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27]
19	Heptachlor	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,9,22] 2) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,27] 3) Soxhlet Extraction, Gas Chromatographic Method ^[10,22] 4) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27]
20	Lead	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^[1,6,15] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14]

3) Digestion...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
21	Lindane	3) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 4) Digestion, Inductively Coupled Plasma Method ^[7,14] 1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,9,22] 2) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,27] 3) Soxhlet Extraction, Gas Chromatographic Method ^[10,22] 4) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27]
22	Mercury	1) Waste Extraction, Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[1,18] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 3) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[19] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]
23	Methoxychlor	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,9,22] 2) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,27] 3) Soxhlet Extraction, Gas Chromatographic Method ^[10,22] 4) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27]


24 Molybdenum...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
24	Molybdenum	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
25	Nickel	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^[1,6,15] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 3) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]
26	Polychlorinated Biphenyls - Aroclor 1016 - Aroclor 1221 - Aroclor 1232 - Aroclor 1242 - Aroclor 1248 - Aroclor 1254 - Aroclor 1260	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[1,9,23] 2) Soxhlet Extraction, Gas Chromatographic Method ^[10,23]
27	Pentachlorophenol	1) Waste Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,25] 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[25]
28	pH	Electrometric Method ^[31,32]
29	Selenium	1) Waste Extraction, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[1,6,20] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 3) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[7,20]


4) Digestion...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
30	Silver	4) Digestion, Inductively Coupled Plasma Method ^[7,14] 1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
31	Thallium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
32	Trichloroethylene	1) Waste Extraction, Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[1,12,26] 2) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[12,26]
33	Vanadium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
34	Zinc	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^[1,6,15] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,14] 3) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]

ดิน จำนวน 124 รายการ

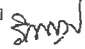ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Acenaphthene	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27] 

2 Acetone...

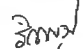ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
2	Acetone	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,26]
3	Aldrin	1) Ultrasonic Extraction, Gas Chromatographic Method ^[11,22] 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,27]
4	Anthracene	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27]
5	Antimony	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[7,16] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
6	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[7,16] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
7	Atrazine	Ultrasonic Extraction, Gas Chromatographic Method ^[11,24]
8	Barium	1) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
9	Benz(a)anthracene	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27]
10	Benzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,26]
11	Benzo(b)fluoranthene	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27]
12	Benzo(k)fluoranthene	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27]
13	Benzoic acid	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,27] 

14 Benzo(a)pyrene...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
14	Benzo(a)pyrene	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,27)
15	Benzo(g,h,i)perylene	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,27)
16	Beryllium	Digestion, Inductively Coupled Plasma Method ^(7,14)
17	Bis(2-chloroethyl)ether	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,27)
18	Bis(2-ethylhexyl)phthalate	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,27)
19	Bromodichloromethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(13,26)
20	Bromoform	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(13,26)
21	Butanol	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(13,26)
22	Butyl benzyl phthalate	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,27)
23	Cadmium	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
24	Carbazole	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,27)
25	Carbon disulfide	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(13,26)
26	Carbon tetrachloride	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(13,26)
27	Chlordane	1) Ultrasonic Extraction, Gas Chromatographic Method ^(11,22) 2) Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,27)

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
28	p-Chloroaniline	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,27)
29	Chlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(13,26)
30	Chlorodibromomethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(13,26)
31	Chloroform	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,26)
32	2-Chlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
33	Chromium	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
34	Chromium (III)	1) Digestion, Flame Atomic Absorption Spectrometric Method; Colorimetric Method; Calculation ^(7,8,15,17) 2) Digestion, Inductively Coupled Plasma Method; Colorimetric Method; Calculation ^(7,8,14,17)
35	Chromium (VI)	Alkaline Digestion, Colorimetric Method ^(8,17)
36	Chrysene	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,27)
37	Cyanide	1) Extraction, Distillation, Titrimetric Method ^(28,29,30) 2) Extraction, Distillation, Colorimetric Method ^(28,29,30)
38	2,4-D	Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽²⁴⁾
39	DDD	1) Ultrasonic Extraction, Gas Chromatographic Method ^(11,22) 2) Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,27)

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
40	DDE	1) Ultrasonic Extraction, Gas Chromatographic Method ^(11,22) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
41	DDT	1) Ultrasonic Extraction, Gas Chromatographic Method ^(11,22) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
42	Dibenz(a,h)anthracene	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,27)
43	Di-n-butyl phthalate	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,27)
44	1,2-Dichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,26)
45	1,3-Dichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,26)
46	1,4-Dichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,26)
47	3,3'-Dichlorobenzidine	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,27)
48	1,1-Dichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,26)
49	1,2-Dichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,26)
50	1,1-Dichloroethylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,26)
51	cis-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,26)
52	trans-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,26)
53	2,4-Dichlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(11,27)

54 1,2-Dichloropropane...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
54	1,2-Dichloropropane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,26)
55	1,3-Dichloropropane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,26)
56	1,3-Dichloropropene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,26)
57	Dieldrin	1) Ultrasonic Extraction, Gas Chromatographic Method ^(11,22) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
58	Diethyl phthalate	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,27)
59	2,4-Dimethylphenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
60	2,4-Dinitrophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
61	2,4-Dinitrotoluene	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,27)
62	2,6-Dinitrotoluene	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,27)
63	Di-n-Octyl phthalate	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,27)
64	Endosulfan	1) Ultrasonic Extraction, Gas Chromatographic Method ^(11,22) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
65	Endrin	1) Ultrasonic Extraction, Gas Chromatographic Method ^(11,22) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
66	Ethylbenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,26)

67 Fluoranthene...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
67	Fluoranthene	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,27)
68	Fluorene	Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,27)
69	Heptachlor	1) Ultrasonic Extraction, Gas Chromatographic Method ^(11,22) 2) Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,27)
70	Heptachlor epoxide	1) Ultrasonic Extraction, Gas Chromatographic Method ^(11,22) 2) Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,27)
71	Hexachlorobenzene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
72	Hexachloro-1,3-butadiene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,26)
73	n-Hexane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(13,26)
74	α-HCH	1) Ultrasonic Extraction, Gas Chromatographic Method ^(11,22) 2) Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,27)
75	β-HCH	1) Ultrasonic Extraction, Gas Chromatographic Method ^(11,22) 2) Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,27)
76	γ-HCH	1) Ultrasonic Extraction, Gas Chromatographic Method ^(11,22) 2) Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,27)
77	Hexachlorocyclopentadiene	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,27)

78 Hexachloroethane...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
78	Hexachloroethane	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,27)
79	Indeno(1,2,3-cd)pyrene	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,27)
80	Isophorone	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,27)
81	Lead	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
82	Manganese	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
83	Mercury	1) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ⁽¹⁹⁾ 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
84	Methanol	Ultrasonic Extraction, Direct Aqueous Injection, Gas Chromatographic Method ^(11,21)
85	Methoxychlor	1) Ultrasonic Extraction, Gas Chromatographic Method ^(11,22) 2) Ultrasonic Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,27)
86	Methyl bromide	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,26)
87	Methylene chloride	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,26)
88	2-Methylphenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
89	2-Methylnaphthalene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(11,27)

90 Methyl tert-butyl ether...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
90	Methyl tert-butyl ether	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,26]
91	Naphthalene	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27]
92	Nickel	1) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
93	Nitrobenzene	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27]
94	N-Nitrosodiphenylamine	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27]
95	N-Nitrosodi-n-propylamine	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27]
96	Polychlorinated Biphenyls - Aroclor 1016 - Aroclor 1221 - Aroclor 1232 - Aroclor 1242 - Aroclor 1248 - Aroclor 1254 - Aroclor 1260	Soxhlet Extraction, Gas Chromatographic Method ^[10,23]
97	Pentachlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[24]
98	Phenanthrene	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27]
99	Phenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,27]
100	Pyrene	Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,27]
101	Selenium	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[7,20]

2) Digestion...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
102	Silver	2) Digestion, Inductively Coupled Plasma Method ^[7,14] 1) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
103	Styrene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,26]
104	1,1,2,2-Tetrachloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,26]
105	Tetrachloroethylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,26]
106	Toluene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,26]
107	TPH (C ₅ -C ₉)	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,26]
108	TPH (C ₈ -C ₁₆)	1) Soxhlet Extraction, Gas Chromatographic Method ^[10,21] 2) Soxhlet Extraction, Gas Chromatographic/Mass spectrometric Method ^[10,26]
109	TPH (C ₁₅ -C ₃₅)	1) Soxhlet Extraction, Gas Chromatographic Method ^[10,21] 2) Soxhlet Extraction, Gas Chromatographic/Mass spectrometric Method ^[10,26]
110	1,2,4-Trichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,26]
111	1,1,1-Trichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,26]
112	1,1,2-Trichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,26]
113	Trichloroethylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,26]

114 2,4,5-Trichlorophenol...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
114	2,4,5-Trichlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,27]
115	2,4,6-Trichlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,27]
116	1,3,5-Trimethylbenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,26]
117	Vanadium	Digestion, Inductively Coupled Plasma Method ^[7,14]
118	Vinyl acetate	Purge and Trap, Gas Chromatographic/Mass spectrometric Method ^[13,26]
119	Vinyl chloride	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,26]
120	m-Xylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,26]
121	o-Xylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,26]
122	p-Xylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,26]
123	Xylene (Total)	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,26]
124	Zinc	1) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 2) Digestion, Inductively Coupled Plasma Method ^[7,14] <i>สิงห์</i>

เอกสารอ้างอิง

- กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2548. เรื่อง การกำจัดสิ่งปฏิกูลหรือวัสดุที่ไม่ใช้แล้ว.ราชกิจจานุเบกษา. 25 มกราคม 2549. เล่มที่ 123 ตอนพิเศษ 11ง.
- กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2549. เรื่อง กำหนดค่าปริมาณเขม่าควันที่เจือปนในอากาศที่ระบายออกจากปล่องของหม้อน้ำโรงสีข้าวที่ใช้กลบเป็นเชื้อเพลิง.ราชกิจจานุเบกษา. 4 ธันวาคม 2549. เล่มที่ 123 ตอนพิเศษ 125ง.
- สมาคมวิศวกรรมสิ่งแวดล้อมแห่งประเทศไทย. คู่มือวิเคราะห์น้ำเสีย. พิมพ์ครั้งที่ 4. กรุงเทพฯ: เรือนแก้วการพิมพ์, 2547.

- APHA, AWWA, WEF. Standard Methods for the Examination of Water and Wastewater. 23rd ed. Washington, DC: APHA, 2017.
- United States Environmental Protection Agency. Standards of Performance for New Stationary Sources. 40 CFR 60. Appendix A, 2023.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. SW-846, 2020.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Acid Digestion of Sediments, Sludges, and Soils. SW-846 Method 3050B, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Alkaline Digestion for Hexavalent Chromium. SW-846 Method 3060A, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Separatory Funnel Liquid-Liquid Extraction. SW-846 Method 3510C, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Soxhlet Extraction. SW-846 Method 3540C, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Ultrasonic Extraction. SW-846 Method 3550C, 2007.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Purge-and-Trap for Aqueous Samples. SW-846 Method 5030C, 2003.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Closed-System Purge-and-Trap and Extraction for Volatile Organics in Soil and Waste Samples. SW-846 Method 5035, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Inductively Coupled Plasma-optical Emission Spectrometry. SW-846 Method 6010D, 2018.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Flame Atomic Absorption Spectrophotometry. SW-846 Method 7000B, 2007.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Antimony and Arsenic (Atomic Absorption, Borohydride Reduction). SW-846 Method 7062, 1994. *สิงห์*

17. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Chromium, Hexavalent (Colorimetric), SW-846 Method 7196A, 1992.

18. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Mercury in Liquid Waste (Manual Cold-Vapor Technique, SW-846 Method 7470A, 1994.

19. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Mercury in Solid or Semisolid Waste (Manual Cold-Vapor Technique, SW-846 Method 7471B, 2007.

20. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Selenium (Atomic Absorption, Borohydride Reduction), SW-846 Method 7742, 1994.

21. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Nonhalogenated Organics Using GC/FID. SW-846 Method 8015D, 2003.


22. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Organochlorine Pesticide by Gas Chromatography. SW-846 Method 8081B, 2007.

23. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Polychlorinated Biphenyls (PCBs) By Gas Chromatography. SW-846 Method 8082A, 2007.

24. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Organophosphorus Compounds by Gas Chromatography. SW-846 Method 8141B, 2007.

25. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Chlorinated Herbicides By GC Using Methylation or Pentafluorobenzoylation Derivatization. SW-846 Method 8151A, 1996.

26. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Volatile Organic Compounds by Gas Chromatography/ Mass Spectrometry (GC/MS). SW-846 Method 8260D, 2018.

27. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. SemiVolatile Organic Compounds by Gas Chromatography/Mass Spectrometry. SW-846 Method 8270E, 2018. 

28. United States...

28. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Total and Amenable Cyanide: Distillation. SW-846 Method 9010C, 2004.

29. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Cyanide Extraction Procedure for Solids and Oils. SW-846 Method 9013A, 2014.

30. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Cyanide in Waters and Extracts Using Titrimetric and Manual Spectrophotometric. SW-846 Method 9014, 2014.

31. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. pH Electrometric Measurement. SW-846 Method 9040C, 2004.

32. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Solid and Waste pH. SW-846 Method 9045D, 2004. 

ภาคผนวก ข

ใบรับรองความสามารถห้องปฏิบัติการและขอข่ายการรับรอง
ห้องปฏิบัติการทดสอบ ตามมาตรฐาน ISO/IEC 17025
จากสำนักงานมาตรฐานอุตสาหกรรม (สมอ.)



แบบ กมช./สมอ.๒
Form NSC/TISI 2

ใบรับรองเลขที่ 24-LB0026
(Certificate No.)

ใบรับรองระบบงาน

(Certificate of Accreditation)

อาศัยอำนาจตามความในพระราชบัญญัติการมาตรฐานแห่งชาติ พ.ศ. ๒๕๕๑
(By Virtue of National Standardization Act B.E. 2551 (2008))

เลขาธิการสำนักงานมาตรฐานผลิตภัณฑ์อุตสาหกรรม
(Secretary-General, Thai Industrial Standards Institute)

ออกใบรับรองฉบับนี้ให้
(Issues this certificate to)

บริษัท ซีคोट จำกัด ฝ่ายห้องปฏิบัติการทดสอบด้านสิ่งแวดล้อม
(Secot Company Limited, Environmental Laboratory Division)

ตั้งอยู่เลขที่
(Address)

๒๓๙ ถนนริมคลองประปา แขวงบางซื่อ เขตบางซื่อ กรุงเทพมหานคร
(239 Rimklongprapa Road, Bangsue, Bangsue, Bangkok)

ได้รับการรับรองความสามารถ
(Certificate of competence)

ตามมาตรฐานเลขที่ มอก. ๑๗๐๒๕ - ๒๕๖๑
(Standard No. TIS 17025-2561 (2018) (ISO/IEC 17025: 2017))

ข้อกำหนดทั่วไปว่าด้วยความสามารถของ ห้องปฏิบัติการทดสอบและห้องปฏิบัติการสอบเทียบ
(General requirements for the competence of testing and calibration laboratories)

หมายเลขการรับรองที่ ทดสอบ ๐๓๙๔
(Accreditation No. Testing 0394)

โดยมีรายละเอียดสาขาและขอบข่ายที่ได้ใบรับรอง แสดงไว้ใน QR CODE และ www.tisi.go.th
(Details of the scheme and scope of the certificate are shown in QR CODE and www.tisi.go.th)

ออกให้ ณ วันที่ ๖ ธันวาคม พ.ศ. ๒๕๖๖
(Issue date : 6 December B.E. 2566 (2023))


(นายวีระศักดิ์ เพ็งหล่ง)
(นายวีระศักดิ์ เพ็งหล่ง)

ผู้อำนวยการสำนักงานคณะกรรมการการมาตรฐานแห่งชาติ
ปฏิบัติราชการแทน
เลขาธิการสำนักงานมาตรฐานผลิตภัณฑ์อุตสาหกรรม



Signed by สำนักงานมาตรฐานผลิตภัณฑ์อุตสาหกรรม (สมอ.)
Thai Industrial Standards Institute (TISI)
Date: 2023-12-06T08:49:54.476+07:00
d68cbe6b

กระทรวงอุตสาหกรรม สำนักงานมาตรฐานผลิตภัณฑ์อุตสาหกรรม
(Ministry of Industry Thailand, Thai Industrial Standards Institute)



รายละเอียดสาขาและขอบข่ายใบรับรองห้องปฏิบัติการ
(Scope of Accreditation for Testing)
ใบรับรองเลขที่ 24-LB0026
(Certification No. 24-LB0026)



ชื่อห้องปฏิบัติการ
(Laboratory Name)

หมายเลขการรับรองที่
(Accreditation No.)

ฉบับที่ 02
(Issue No.02)

สถานภาพห้องปฏิบัติการ
(Laboratory status)

บริษัท ซีคोट จำกัด ฝ่ายห้องปฏิบัติการทดสอบด้านสิ่งแวดล้อม
(Secot Company Limited, Environmental Laboratory Division)

ทดสอบ 0394
(Testing 0394)

ออกให้ตั้งแต่วันที่ 30 ตุลาคม พ.ศ. 2566
(Valid from) (30 October B.E.2566 (2023))

☒ ถาวร (Permanent) ☐ นอกสถานที่ (Site) ☐ชั่วคราว (Temporary)

ถึงวันที่ 8 กันยายน พ.ศ. 2571
(Until) (8 September B.E.2571 (2028))

☐เคลื่อนที่ (Mobile) ☐หลายสถานที่ (Multisite)

สาขาการทดสอบ (Field of Testing)	รายการทดสอบ (Parameter)	วิธีทดสอบ (Test Method)
สาขาสสิ่งแวดล้อม (environmental field) 1. น้ำและน้ำเสีย (water and wastewater)	- โลหะหนัก (heavy metals) • สารหนู (Arsenic, As) 0.000 5 mg/L ถึง 0.090 0 mg/L • สารหนู (Arsenic, As) 0.05 mg/L ถึง 4.50 mg/L • แบเรียม (Barium, Ba) 0.02 mg/L ถึง 4.50 mg/L • แคดเมียม (Cadmium, Cd) 0.01 mg/L ถึง 4.50 mg/L • โครเมียม (Chromium, Cr) 0.01 mg/L ถึง 4.50 mg/L	- Standard Methods for the Examination of Water and Wastewater, APHA, AWWA, WEF, 23 rd edition, 2017, Part 3030 F and Part 3114 C - Standard Methods for the Examination of Water and Wastewater, APHA, AWWA, WEF, 23 rd edition, 2017, Part 3030 E and Part 3120 B

กระทรวงอุตสาหกรรมสำนักงานมาตรฐานผลิตภัณฑ์อุตสาหกรรม
(Ministry of Industry, Thai Industrial Standards Institute)

หน้า 1/9

รายละเอียดสาขาและขอบข่ายใบรับรองห้องปฏิบัติการ

(Scope of Accreditation for Testing)

ใบรับรองเลขที่ 24-LB0026

(Certification No. 24-LB0026)



ฉบับที่ 02

(Issue No.02)

ออกให้ตั้งแต่วันที่ 30 ตุลาคม พ.ศ. 2566

(Valid from)

(30 October B.E.2566 (2023))

ถึงวันที่ 8 กันยายน พ.ศ. 2571

(Until) (8 September B.E.2571 (2028))

สถานภาพห้องปฏิบัติการ
(Laboratory status)

☒ ถาวร
(Permanent)

☐ นอกสถานที่
(Site)

☐ชั่วคราว
(Temporary)

☐เคลื่อนที่
(Mobile)

☐หลายสถานที่
(Multisite)

สาขาการทดสอบ (Field of Testing)	รายการทดสอบ (Parameter)	วิธีทดสอบ (Test Method)
สาขาสิ่งแวดล้อม (environmental field)		
1. น้ำและน้ำเสีย (ต่อ) (water and wastewater) (cont.)	<ul style="list-style-type: none"> - โลหะหนัก (heavy metals) • ทองแดง (Copper, Cu) 0.02 mg/L ถึง 4.50 mg/L • เหล็ก (Iron, Fe) 0.05 mg/L ถึง 9.00 mg/L • ตะกั่ว (Lead, Pb) 0.03 mg/L ถึง 4.50 mg/L • แมงกานีส (Manganese, Mn) 0.01 mg/L ถึง 9.00 mg/L • นิกเกิล (Nickel, Ni) 0.01 mg/L ถึง 4.50 mg/L • สังกะสี (Zinc, Zn) 0.02 mg/L ถึง 9.00 mg/L 	<ul style="list-style-type: none"> - Standard Methods for the Examination of Water and Wastewater, APHA, AWWA, WEF, 23rd edition, 2017, Part 3030 E and Part 3120 B

กระทรวงอุตสาหกรรมสำนักงานมาตรฐานผลิตภัณฑ์อุตสาหกรรม

(Ministry of Industry, Thai Industrial Standards Institute)

หน้าที่ 2/9

รายละเอียดสาขาและขอบข่ายใบรับรองห้องปฏิบัติการ

(Scope of Accreditation for Testing)

ใบรับรองเลขที่ 24-LB0026

(Certification No. 24-LB0026)



ฉบับที่ 02

(Issue No.02)

ออกให้ตั้งแต่วันที่ 30 ตุลาคม พ.ศ. 2566

(Valid from)

(30 October B.E.2566 (2023))

ถึงวันที่ 8 กันยายน พ.ศ. 2571

(Until) (8 September B.E.2571 (2028))

สถานภาพห้องปฏิบัติการ
(Laboratory status)

☒ ถาวร
(Permanent)

☐นอกสถานที่
(Site)

☐ชั่วคราว
(Temporary)

☐เคลื่อนที่
(Mobile)

☐หลายสถานที่
(Multisite)

สาขาการทดสอบ (Field of Testing)	รายการทดสอบ (Parameter)	วิธีทดสอบ (Test Method)
สาขาสิ่งแวดล้อม (environmental field)		
1. น้ำและน้ำเสีย (ต่อ) (water and wastewater) (cont.)	<ul style="list-style-type: none"> - ซีโอดี (Chemical oxygen demand, COD) 100 mg/L ถึง 4 000 mg/L 	<ul style="list-style-type: none"> - Standard Methods for the Examination of Water and Wastewater, APHA, AWWA, WEF, 23rd edition, 2017, Part 5220 D
2. บริเวณทำงาน (workplace)	<ul style="list-style-type: none"> - ฝุ่นละอองรวม (Total dust) 0.10 mg/filter ถึง 2.00 mg/filter - ฝุ่นละอองขนาดเล็ก (Respirable dust) 0.10 mg/filter ถึง 2.00 mg/filter 	<ul style="list-style-type: none"> - NIOSH Manual of Analytical Methods (NMAM), method 0500, 4th edition, 15th August 1994 (Exclude Sampling) - NIOSH Manual of Analytical Methods (NMAM), method 0600, 4th edition, 15th January 1998 (Exclude Sampling)

กระทรวงอุตสาหกรรมสำนักงานมาตรฐานผลิตภัณฑ์อุตสาหกรรม

(Ministry of Industry, Thai Industrial Standards Institute)

หน้าที่ 3/9

รายละเอียดสาขาและขอบข่ายใบรับรองห้องปฏิบัติการ

(Scope of Accreditation for Testing)

ใบรับรองเลขที่ 24-LB0026

(Certification No. 24-LB0026)



ฉบับที่ 02
(Issue No.02)

ออกให้ตั้งแต่วันที่ 30 ตุลาคม พ.ศ. 2566
(Valid from) (30 October B.E.2566 (2023))

ถึงวันที่ 8 กันยายน พ.ศ. 2571
(Until) (8 September B.E.2571 (2028))

สถานภาพห้องปฏิบัติการ
(Laboratory status)

☒ ถาวร
(Permanent)

☐ นอกสถานที่
(Site)

☐ชั่วคราว
(Temporary)

☐เคลื่อนที่
(Mobile)

☐หลายสถานที่
(Multisite)

สาขาการทดสอบ (Field of Testing)	รายการทดสอบ (Parameter)	วิธีทดสอบ (Test Method)
<p>สาขาส่งแวดล้อม (environmental field)</p> <p>2. บริเวณทำงาน (ต่อ) (workplace) (cont.)</p>	<ul style="list-style-type: none"> เบนซีน (Benzene) 1.10 µg/tube ถึง 420 µg/tube โทลูอีน (Toluene) 1.10 µg/tube ถึง 420 µg/tube โทโครไซลีน (Total xylenes) 2.20 µg/tube ถึง 840 µg/tube เมตา, พารา-ไซลีน (m, p- Xylene) 1.10 µg/tube ถึง 420 µg/tube ออร์โธ-ไซลีน (o- Xylene) 1.10 µg/tube ถึง 420 µg/tube 	<ul style="list-style-type: none"> NIOSH Manual of Analytical Methods (NMAM) , method 1501, 4th edition , 15th March 2003 (Exclude Sampling)
<p>3. ปล่องระบายอากาศ (stack)</p>	<ul style="list-style-type: none"> ซัลเฟอร์ไดออกไซด์ (Sulfur dioxide) 1.00 mg/L ถึง 16 000 mg/L (solution) 	<ul style="list-style-type: none"> US.EPA , Code of Federal Regulations , 40 CFR 60 appendix A , method 6 , July 2019 (Exclude Sampling)

รายละเอียดสาขาและขอบข่ายใบรับรองห้องปฏิบัติการ

(Scope of Accreditation for Testing)

ใบรับรองเลขที่ 24-LB0026

(Certification No. 24-LB0026)



ฉบับที่ 02
(Issue No.02)

ออกให้ตั้งแต่วันที่ 30 ตุลาคม พ.ศ. 2566
(Valid from) (30 October B.E.2566 (2023))

ถึงวันที่ 8 กันยายน พ.ศ. 2571
(Until) (8 September B.E.2571 (2028))

สถานภาพห้องปฏิบัติการ
(Laboratory status)

☒ ถาวร
(Permanent)

☐นอกสถานที่
(Site)

☐ชั่วคราว
(Temporary)

☐เคลื่อนที่
(Mobile)

☐หลายสถานที่
(Multisite)

สาขาการทดสอบ (Field of Testing)	รายการทดสอบ (Parameter)	วิธีทดสอบ (Test Method)
<p>สาขาส่งแวดล้อม (environmental field)</p> <p>3. ปล่องระบายอากาศ (ต่อ) (stack) (cont.)</p>	<ul style="list-style-type: none"> ไฮโดรเจนฟลูออไรด์ (Hydrogen fluoride) 5 µg/sample ถึง 400 µg/sample ไฮโดรเจนคลอไรด์ (Hydrogen chloride) 5 µg/sample ถึง 400 µg/sample 	<ul style="list-style-type: none"> WI-7.2-1-22 based on US.EPA , Code of Federal Regulations , 40 CFR 60 appendix A , method 26 , 2019 (Exclude Sampling)

รายละเอียดสาขาและขอบข่ายใบรับรองห้องปฏิบัติการ

(Scope of Accreditation for Testing)

ใบรับรองเลขที่ 24-LB0026

(Certification No. 24-LB0026)



ฉบับที่ 02

(Issue No.02)

ออกให้ตั้งแต่วันที่ 30 ตุลาคม พ.ศ. 2566

(Valid from)

(30 October B.E.2566 (2023))

ถึงวันที่ 8 กันยายน พ.ศ. 2571

(Until) (8 September B.E.2571 (2028))

สถานภาพห้องปฏิบัติการ
(Laboratory status)

☒ถาวร
(Permanent)

☒นอกสถานที่
(Site)

☐ชั่วคราว
(Temporary)

☐เคลื่อนที่
(Mobile)

☐หลายสถานที่
(Multisite)

สาขาการทดสอบ (Field of Testing)	รายการทดสอบ (Parameter)	วิธีทดสอบ (Test Method)
<p>สาขาส่งแวดล้อม (environmental field)</p> <p>4. บรรยากาศทั่วไป (ambient air)</p>	<p>- สารอินทรีย์ระเหยง่าย (Volatile organic compounds, VOCs)</p> <ul style="list-style-type: none"> คลอโรอีเทน (Chloroethene) 0.05 $\mu\text{g}/\text{m}^3$ ถึง 51.00 $\mu\text{g}/\text{m}^3$ (0.02 ppbv ถึง 20.00 ppbv) 1,3-บิวทาไดเ็น (1,3-butadiene) 0.04 $\mu\text{g}/\text{m}^3$ ถึง 44.00 $\mu\text{g}/\text{m}^3$ (0.02 ppbv ถึง 20.00 ppbv) โบรมอมีเทน (Bromomethane) 0.08 $\mu\text{g}/\text{m}^3$ ถึง 77.00 $\mu\text{g}/\text{m}^3$ (0.02 ppbv ถึง 20.00 ppbv) อะคลอลีน (Acrolein) 0.05 $\mu\text{g}/\text{m}^3$ ถึง 45.00 $\mu\text{g}/\text{m}^3$ (0.02 ppbv ถึง 20.00 ppbv) 	<p>- WI-7.2-1-24 based on US EPA , Compendium Method TO-15 , EPA/625/R-96/010b, Second edition, January 1999</p>

กระทรวงอุตสาหกรรมสำนักงานมาตรฐานผลิตภัณฑ์อุตสาหกรรม

(Ministry of Industry, Thai Industrial Standards Institute)

หน้า 6/9

รายละเอียดสาขาและขอบข่ายใบรับรองห้องปฏิบัติการ

(Scope of Accreditation for Testing)

ใบรับรองเลขที่ 24-LB0026

(Certification No. 24-LB0026)



ฉบับที่ 02

(Issue No.02)

ออกให้ตั้งแต่วันที่ 30 ตุลาคม พ.ศ. 2566

(Valid from)

(30 October B.E.2566 (2023))

ถึงวันที่ 8 กันยายน พ.ศ. 2571

(Until) (8 September B.E.2571 (2028))

สถานภาพห้องปฏิบัติการ
(Laboratory status)

☒ถาวร
(Permanent)

☒นอกสถานที่
(Site)

☐ชั่วคราว
(Temporary)

☐เคลื่อนที่
(Mobile)

☐หลายสถานที่
(Multisite)

สาขาการทดสอบ (Field of Testing)	รายการทดสอบ (Parameter)	วิธีทดสอบ (Test Method)
<p>สาขาส่งแวดล้อม (environmental field)</p> <p>4. บรรยากาศทั่วไป (ต่อ) (ambient air) (cont.)</p>	<p>- สารอินทรีย์ระเหยง่าย (Volatile organic compounds, VOCs)</p> <ul style="list-style-type: none"> อะครีโลไนไทรล์ (Acrylonitrile) 0.04 $\mu\text{g}/\text{m}^3$ ถึง 43.00 $\mu\text{g}/\text{m}^3$ (0.02 ppbv ถึง 20.00 ppbv) ไดคลอโรมีเทน (Dichloromethane) 0.14 $\mu\text{g}/\text{m}^3$ to 69.00 $\mu\text{g}/\text{m}^3$ (0.04 ppbv ถึง 20.00 ppbv) คาร์บอนไดซัลไฟด์ (Carbon disulfide) 0.06 $\mu\text{g}/\text{m}^3$ ถึง 62.00 $\mu\text{g}/\text{m}^3$ (0.02 ppbv ถึง 20.00 ppbv) ไตรคลอโรมีเทน (Trichloromethane) 0.20 $\mu\text{g}/\text{m}^3$ ถึง 97.00 $\mu\text{g}/\text{m}^3$ (0.04 ppbv ถึง 20.00 ppbv) 1,2-ไดคลอโรอีเทน (1,2-dichloroethane) 0.08 $\mu\text{g}/\text{m}^3$ ถึง 80.00 $\mu\text{g}/\text{m}^3$ (0.02 ppbv ถึง 20.00 ppbv) 	<p>- WI-7.2-1-24 based on US EPA , Compendium Method TO-15 , EPA/625/R-96/010b, Second edition, January 1999</p>

กระทรวงอุตสาหกรรมสำนักงานมาตรฐานผลิตภัณฑ์อุตสาหกรรม

(Ministry of Industry, Thai Industrial Standards Institute)

หน้า 7/9

รายละเอียดสาขาและขอบข่ายใบรับรองห้องปฏิบัติการ

(Scope of Accreditation for Testing)

ใบรับรองเลขที่ 24-LB0026

(Certification No. 24-LB0026)



ฉบับที่ 02
(Issue No.02)

ออกให้ตั้งแต่วันที่ 30 ตุลาคม พ.ศ. 2566
(Valid from) (30 October B.E.2566 (2023))

ถึงวันที่ 8 กันยายน พ.ศ. 2571
(Until) (8 September B.E.2571 (2028))

สถานภาพห้องปฏิบัติการ
(Laboratory status)

☒ ถาวร
(Permanent)

☒ นอกสถานที่
(Site)

☐ชั่วคราว
(Temporary)

☐เคลื่อนที่
(Mobile)

☐หลายสถานที่
(Multisite)

สาขาการทดสอบ (Field of Testing)	รายการทดสอบ (Parameter)	วิธีทดสอบ (Test Method)
<p>สาขาส่งแวดล้อม (environmental field)</p> <p>4. บรรยากาศทั่วไป (ต่อ) (ambient air) (cont.)</p>	<p>- สารอินทรีย์ระเหยง่าย (Volatile organic compounds, VOCs)</p> <ul style="list-style-type: none"> เบนซีน (Benzene) 0.06 $\mu\text{g}/\text{m}^3$ ถึง 63.00 $\mu\text{g}/\text{m}^3$ (0.02 ppbv ถึง 20.00 ppbv) คาร์บอนเตตระคลอไรด์ (Carbon tetrachloride) 0.25 $\mu\text{g}/\text{m}^3$ ถึง 125 $\mu\text{g}/\text{m}^3$ (0.04 ppbv ถึง 20.00 ppbv) ไตรคลอโรเอทิลีน (Trichloroethylene) 0.21 $\mu\text{g}/\text{m}^3$ ถึง 107 $\mu\text{g}/\text{m}^3$ (0.04 ppbv ถึง 20.00 ppbv) 1,2-ไดคลอโรโพรเพน (1,2-dichloropropane) 0.18 $\mu\text{g}/\text{m}^3$ ถึง 92.00 $\mu\text{g}/\text{m}^3$ (0.04 ppbv ถึง 20.00 ppbv) เตตระคลอโรเอทิลีน (Tetrachloroethylene) 0.27 $\mu\text{g}/\text{m}^3$ ถึง 135 $\mu\text{g}/\text{m}^3$ (0.04 ppbv ถึง 20.00 ppbv) 	<p>- WI-7.2-1-24 based on US EPA , Compendium Method TO-15 , EPA/625/R-96/010b, Second edition, January 1999</p>

รายละเอียดสาขาและขอบข่ายใบรับรองห้องปฏิบัติการ

(Scope of Accreditation for Testing)

ใบรับรองเลขที่ 24-LB0026

(Certification No. 24-LB0026)



ฉบับที่ 02
(Issue No.02)

ออกให้ตั้งแต่วันที่ 30 ตุลาคม พ.ศ. 2566
(Valid from) (30 October B.E.2566 (2023))

ถึงวันที่ 8 กันยายน พ.ศ. 2571
(Until) (8 September B.E.2571 (2028))

สถานภาพห้องปฏิบัติการ
(Laboratory status)

☒ ถาวร
(Permanent)

☒ นอกสถานที่
(Site)

☐ชั่วคราว
(Temporary)

☐เคลื่อนที่
(Mobile)

☐หลายสถานที่
(Multisite)

สาขาการทดสอบ (Field of Testing)	รายการทดสอบ (Parameter)	วิธีทดสอบ (Test Method)
<p>สาขาส่งแวดล้อม (environmental field)</p> <p>4. บรรยากาศทั่วไป (ต่อ) (ambient air) (cont.)</p>	<p>- สารอินทรีย์ระเหยง่าย (Volatile organic compounds ,VOCs)</p> <ul style="list-style-type: none"> 1,2-ไดโบรมโอเอเทน (1,2-dibromoethane) 0.31 $\mu\text{g}/\text{m}^3$ ถึง 153 $\mu\text{g}/\text{m}^3$ (0.04 ppbv ถึง 20.00 ppbv) 1,1,2,2-เตตระคลอโรอีเทน (1,1,2,2-tetrachloroethane) 0.69 $\mu\text{g}/\text{m}^3$ ถึง 137 $\mu\text{g}/\text{m}^3$ (0.10 ppbv ถึง 20.00 ppbv) เบนซิลคลอไรด์ (Benzyl chloride) 0.52 $\mu\text{g}/\text{m}^3$ ถึง 103 $\mu\text{g}/\text{m}^3$ (0.10 ppbv ถึง 20.00 ppbv) 1,4-ไดคลอโรเบนซีน (1,4-dichlorobenzene) 0.24 $\mu\text{g}/\text{m}^3$ ถึง 120 $\mu\text{g}/\text{m}^3$ (0.04 ppbv ถึง 20.00 ppbv) 	<p>- WI-7.2-1-24 based on US EPA , Compendium Method TO-15 , EPA/625/R-96/010b, Second edition, January 1999</p>

ภาคผนวก ซ

ใบอนุญาตเป็นนิติบุคคลผู้ให้บริการตรวจวัดและวิเคราะห์การทำงาน
จากกรมสวัสดิการและคุ้มครองแรงงาน



แบบ กภ.บุญ
นิติบุคคล

กรมสวัสดิการและคุ้มครองแรงงาน
ใบอนุญาต
เป็นนิติบุคคลผู้ให้บริการตรวจวัดและวิเคราะห์สภาวะการทำงานเกี่ยวกับระดับเสียง

ใบอนุญาตเลขที่ ๑๔๑๓-๐๓-๒๕๖๕-๐๐๔๘

อนุญาตให้.....บริษัท ชีคอต จั๊กกิต.....

เลขทะเบียนนิติบุคคล.....๐๑๐๕๕๓๖๐๐๐๙๗๖.....

ตั้งอยู่ เลขที่ ๒๓๙ ถนนริมคลองประปา แขวงนางชื่อ เขตนางชื่อ กรุงเทพมหานคร.....

เป็นนิติบุคคลผู้ให้บริการด้านความปลอดภัย อาชีวอนามัย และสภาพแวดล้อมในการทำงาน ตามกฎกระทรวงกำหนดมาตรฐานในการบริหาร จัดการ และดำเนินการด้านความปลอดภัย อาชีวอนามัย และสภาพแวดล้อมในการทำงาน เกี่ยวกับความร้อน แสงสว่าง และเสียง พ.ศ. ๒๕๕๙ ในการตรวจวัดและวิเคราะห์สภาวะการทำงานเกี่ยวกับระดับเสียง ประกอบกับกฎกระทรวงการขึ้นทะเบียนและการอนุญาตให้บริการเพื่อส่งเสริม ความปลอดภัย อาชีวอนามัย และสภาพแวดล้อมในการทำงาน พ.ศ. ๒๕๖๔ แห่งพระราชบัญญัติความปลอดภัย อาชีวอนามัย และสภาพแวดล้อมในการทำงาน พ.ศ. ๒๕๕๔ โดยมีบุคลากร จำนวน ๕ ราย ดังรายชื่อแนบท้ายใบอนุญาตนี้

ทั้งนี้ ตั้งแต่วันที่ ๑๗ มิถุนายน พ.ศ. ๒๕๖๕ ถึงวันที่ ๑๖ มิถุนายน พ.ศ. ๒๕๖๘

ให้ไว้ ณ วันที่ ๑๗ มิถุนายน พ.ศ. ๒๕๖๕

(นายสมพจน์ กวางแก้ว)
รองอธิบดี ปฏิบัติราชการแทน
อธิบดีกรมสวัสดิการและคุ้มครองแรงงาน

เลขทะเบียนควบคุม

๗-๑๑-๐๔๐๓-๐๔๘-๐๑-๖๕

(ลงนาม)

(นายทะเบียน)

(นายศักดิ์ศิลป์ ตูลาธร)

ตำแหน่ง ผู้อำนวยการกองความปลอดภัยแรงงาน

รายชื่อบุคลากรแนบท้ายใบอนุญาต
เป็นนิติบุคคลผู้ให้บริการตรวจวัดและวิเคราะห์สภาพการทำงานเกี่ยวกับระดับเสียง
ของบริษัท ชีคอท จำกัด

ใบอนุญาตเลขที่ ๐๔๐๓-๐๓-๒๕๖๕-๐๐๔๘

- | | |
|-------------------|----------------|
| ๑. นางสาวสุนันทา | ศิริพัฒน์นันท์ |
| ๒. นางสาวกนิษฐา | เจริญเชื้อ |
| ๓. นางสาวปัทมวรรณ | สุวรรณวิโรจน์ |
| ๔. นางสาวอลิษา | คณิธรานนท์ |
| ๕. นางสาวชนิตา | หล้าสาย |

ทั้งนี้ ตั้งแต่วันที่ ๑๗ มิถุนายน พ.ศ. ๒๕๖๕ ถึงวันที่ ๑๖ มิถุนายน พ.ศ. ๒๕๖๘

ให้ไว้ ณ วันที่ ๑๗ มิถุนายน พ.ศ. ๒๕๖๕



(นายสมพนธ์ กวางแก้ว)

รองอธิบดี ปฏิบัติราชการแทน
อธิบดีกรมสวัสดิการและคุ้มครองแรงงาน


รายชื่อบุคลากร (เพิ่มเติม)
แนบท้ายใบอนุญาตเป็นนิติบุคคลผู้ให้บริการตรวจวัดและวิเคราะห์สภาพการทำงานเกี่ยวกับระดับเสียง
ของบริษัท ชีคอท จำกัด

ใบอนุญาตเลขที่ ๐๔๐๓-๐๓-๒๕๖๕-๐๐๔๘

- | | |
|-------------------|-------------|
| ๑. นางสาวศลิษา | อินริย์ |
| ๒. นางสาวมาริยาณี | อาแว |
| ๓. นางสาววิระยา | ปัจฉิมบุรณ์ |

ทั้งนี้ ตั้งแต่วันที่ ๑๓ มกราคม พ.ศ. ๒๕๖๖ ถึงวันที่ ๑๖ มิถุนายน พ.ศ. ๒๕๖๘

ให้ไว้ ณ วันที่ ๑๓ มกราคม พ.ศ. ๒๕๖๖

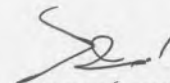

(นายสมพนธ์ กวางแก้ว)

รองอธิบดี ปฏิบัติราชการแทน
อธิบดีกรมสวัสดิการและคุ้มครองแรงงาน



แบบ ก.ก.ม.อ.
อ.ค.อ.

กรมสวัสดิการและคุ้มครองแรงงาน

ใบอนุญาต

เป็นนิติบุคคลผู้ให้บริการตรวจวัดระดับความเข้มข้นของสารเคมีอันตราย
ในบรรยากาศของสถานที่ทำงานและสถานที่เก็บรักษาสารเคมีอันตราย

ใบอนุญาตเลขที่ ๐๒๐๑-๑๓-๒๕๖๕-๐๐๔๙

อนุญาตให้ บริษัท ซีคอน จำกัด

เลขทะเบียนนิติบุคคล ๐๑๐๕๕๓๖๐๐๐๙๗๖

ตั้งอยู่ เลขที่ ๒๓๙ ถนนริมคลองประปา แขวงบางซื่อ เขตบางซื่อ กรุงเทพมหานคร
เป็นนิติบุคคลผู้ให้บริการด้านความปลอดภัย อาชีวอนามัย และสภาพแวดล้อมในการทำงาน ตามกฎกระทรวง
กำหนดมาตรฐานในการบริหาร จัดการ และดำเนินการด้านความปลอดภัย อาชีวอนามัย และสภาพแวดล้อม
ในการทำงานเกี่ยวกับสารเคมีอันตราย พ.ศ. ๒๕๕๖ ในการเป็นผู้ให้บริการตรวจวัดระดับความเข้มข้น
ของสารเคมีอันตรายในบรรยากาศของสถานที่ทำงานและสถานที่เก็บรักษาสารเคมีอันตราย ประกอบกับ
กฎกระทรวงการขึ้นทะเบียนและการอนุญาตให้บริการเพื่อส่งเสริมความปลอดภัย อาชีวอนามัย และสภาพแวดล้อม
ในการทำงาน พ.ศ. ๒๕๖๔ แห่งพระราชบัญญัติความปลอดภัย อาชีวอนามัย และสภาพแวดล้อมในการทำงาน
พ.ศ. ๒๕๕๔ โดยมีบุคลากร จำนวน ๑๔ ราย ดังรายชื่อแนบท้ายใบอนุญาตนี้

ทั้งนี้ ตั้งแต่วันที่ ๑๔ มิถุนายน พ.ศ. ๒๕๖๕ ถึงวันที่ ๑๓ มิถุนายน พ.ศ. ๒๕๖๘

ให้ไว้ ณ วันที่ ๑๔ มิถุนายน พ.ศ. ๒๕๖๕

(นายสมพงษ์ กวางแก้ว)

รองอธิบดี ปฏิบัติราชการแทน

อธิบดีกรมสวัสดิการและคุ้มครองแรงงาน

เลขทะเบียนควบคุม

ข-๑๑-๐๒๐๑-๐๔๙-๐๓-๖๕

(ลงนาม)

(นายทะเบียน)

(นายศักดิ์ศิลป์ ตูลาธร)


ผู้อำนวยการกองความปลอดภัยแรงงาน

รายชื่อบุคลากรแนบท้ายใบอนุญาต
เป็นนิติบุคคลผู้ให้บริการตรวจวัดระดับความเข้มข้นของสารเคมีอันตรายในบรรยากาศของสถานที่ทำงาน
และสถานที่เก็บรักษาสารเคมีอันตราย
ของบริษัท ซีคอต จำกัด
ใบอนุญาตเลขที่ ๐๒๐๑-๐๓-๒๕๖๕-๐๐๔๔

๑. นายชิตพล	สมประสงค์
๒. นายอนิวัฒน์	พิมพ์นา
๓. นายศิวะนนท์	กุลวงษ์
๔. นายวัชรกานต์	ประมาคะเต
๕. นายธนโชติ	ช่างลือ
๖. นายกิตติพงษ์	ละเกิงสุข
๗. นายจิรวีวัฒน์	โคตรคำหาญ
๘. นายศุภกิจ	ติ่มภูกา
๙. นางสาวธัญลักษณ์	โยธา
๑๐. นางสาวทิพย์สุดา	วรรณการ
๑๑. นางสาวสายธาร	ภูเขียว
๑๒. นายภาคภูมิ	แทนไทย
๑๓. นายธนาวุฒิ	ถ้วนแสง
๑๔. นายรัตนชัย	ชอบทำกิจ

ทั้งนี้ ตั้งแต่วันที่ ๑๔ มิถุนายน พ.ศ. ๒๕๖๕ ถึงวันที่ ๑๓ มิถุนายน พ.ศ. ๒๕๖๘

ให้ไว้ ณ วันที่ ๑๔ มิถุนายน พ.ศ. ๒๕๖๕


(นายสมพงษ์ กวางแก้ว)
รองอธิบดี ปฏิบัติราชการแทน
อธิบดีกรมสวัสดิการและคุ้มครองแรงงาน



แบบ กภ.บอญ
ฉศ.บุคคล

กรมสวัสดิการและคุ้มครองแรงงาน
ใบอนุญาต

เป็นนิติบุคคลผู้ให้บริการวิเคราะห์ระดับความเข้มข้นของสารเคมีอันตราย
ในบรรยากาศของสถานที่ทำงานและสถานที่เก็บรักษาสารเคมีอันตราย


ใบอนุญาตเลขที่ ๐๒๐๑-๐๓-๒๕๖๕-๐๐๓๔

อนุญาตให้ บริษัท ซีคอต จำกัด

เลขทะเบียนนิติบุคคล ๐๑๐๕๕๓๖๐๐๐๗๖
ตั้งอยู่ เลขที่ ๒๓๙ ถนนริมคลองประปา แขวงบางซื่อ เขตบางซื่อ กรุงเทพมหานคร
เป็นนิติบุคคลผู้ให้บริการด้านความปลอดภัย อาชีวอนามัย และสภาพแวดล้อมในการทำงาน ตามกฎกระทรวง
กำหนดมาตรฐานในการบริหาร จัดการ และดำเนินการด้านความปลอดภัย อาชีวอนามัย และสภาพแวดล้อม
ในการทำงานเกี่ยวกับสารเคมีอันตราย พ.ศ. ๒๕๕๖ ในการเป็นผู้ให้บริการวิเคราะห์ระดับความเข้มข้น
ของสารเคมีอันตรายในบรรยากาศของสถานที่ทำงานและสถานที่เก็บรักษาสารเคมีอันตราย ประกอบกับ
กฎกระทรวงการขึ้นทะเบียนและการอนุญาตให้บริการเพื่อส่งเสริมความปลอดภัย อาชีวอนามัย และสภาพแวดล้อม
ในการทำงาน พ.ศ. ๒๕๖๔ แห่งพระราชบัญญัติความปลอดภัย อาชีวอนามัย และสภาพแวดล้อมในการทำงาน
พ.ศ. ๒๕๕๔ โดยมีบุคลากร จำนวน ๑๔ ราย ดังรายชื่อแนบท้ายใบอนุญาตนี้

ทั้งนี้ ตั้งแต่วันที่ ๑๔ มิถุนายน พ.ศ. ๒๕๖๕ ถึงวันที่ ๑๓ มิถุนายน พ.ศ. ๒๕๖๘

ให้ไว้ ณ วันที่ ๑๔ มิถุนายน พ.ศ. ๒๕๖๕


(นายสมพงษ์ กวางแก้ว)
รองอธิบดี ปฏิบัติราชการแทน
อธิบดีกรมสวัสดิการและคุ้มครองแรงงาน

รายชื่อบุคลากรแนบท้ายใบอนุญาต
เป็นนิติบุคคลผู้ให้บริการวิเคราะห์ระดับความเข้มข้นของสารเคมีอันตรายในบรรยากาศของสถานที่ทำงาน
และสถานที่เก็บรักษาสารเคมีอันตราย
ของบริษัท ซีคอท จำกัด
ใบอนุญาตเลขที่ ๐๒๐๒-๐๓-๒๖๕-๐๐๓๔

- | | |
|-----------------------|----------------|
| ๑. นางสาวนริสา | ภูวสรเพ็ชญ์ |
| ๒. นางอารยา | ทิทธิรักษ์ |
| ๓. นางสาวศิริวรรณ | นิมสง่า |
| ๔. นางสาวสุชาทิพย์ | เทียนเตี้ย |
| ๕. นางสาวพรนภา | บุศรธรรม |
| ๖. นางสาวชาวีณี | อาจปลิว |
| ๗. นางสาวกฤษณา | จันทุม |
| ๘. นางสาวพัชรา | สมานอินท์ |
| ๙. นางสาวณิสตา | ก๊วยออน |
| ๑๐. นางสาวศศิภา | ใจดี |
| ๑๑. นางสาวจุฑารัตน์ | แจ่มเรือน |
| ๑๒. นางสาวณัฐศิริ | เลิศธีรพัฒน์ |
| ๑๓. นางสาวสัญญาลักษณ์ | อินทรประสิทธิ์ |
| ๑๔. นางสาวสุดาพร | สุนทร |

ทั้งนี้ ตั้งแต่วันที่ ๑๕ มิถุนายน พ.ศ. ๒๕๖๕ ถึงวันที่ ๑๓ มิถุนายน พ.ศ. ๒๕๖๘

ให้ไว้ ณ วันที่ ๑๕ มิถุนายน พ.ศ. ๒๕๖๕



(นายสมพงษ์ ภาวนานันท์)

รองอธิบดี ปฏิบัติราชการแทน

อธิบดีกรมสวัสดิการและคุ้มครองแรงงาน

เลขทะเบียนควบคุม
ข-๑๑-๐๒๐๒-๐๓๔-๐๑-๖๕

(ลงนาม).....(นายทะเบียน)

(นายศักดิ์ศิลป์ ตูลารธร)

ผู้อำนวยการกองความปลอดภัยแรงงาน