

ภาคผนวก ง

เอกสารสอบเทียบเครื่องมือ

สีเหลือง=USER ผู้ร้องขอทำออก

List of Instrument Certificates for Environmental Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*
1	Atomic Absorption Spectrometer	CADMIUM COPPER LEAD MANGANESE MERCURY NICKEL TOTAL IRON ZINC	Agilent Technologies	AA240FS / MY13160001	Thailand Institute of Scientific and Technological Research (TISTR)	MTC ACL No 387 66	2 Feb 23	1 Feb 24
2	Atomic Absorption Spectrometer	ARSENIC SELENIUM	Perkin Elmer	PinAAcle 900F / PFB20031902	Perkin Elmer Co.,Ltd.	WO-02273773	26 Jun 23	25 Jun 24
3	Analytical Balance	FAT OIL AND GREASE	Mettler Toledo	AB204-S/FACT / 1129361010	National Food Institute,Ministry of Industry, Thailand	2303074-001-01	27 May 23	25 May 24
4	Analytical Balance	TOTAL DISSOLVED SOLIDS TOTAL SUSPENDED SOLIDS	Mettler Toledo	XSR205DU / C009071872	Technology Promotion Association (Thailand~Japan)	23MM112	26 Apr 23	25 Apr 24
5	DO Meter	DO	Horiba	LAQUA-DO210 / HE1M0020	Technology Promotion Association (Thailand~Japan)	23TW90	29 Mar 23	28 Mar 24
6	DO Meter	BIOCHEMICAL OXYGEN DEMAND	YSI	5100 / 11B 101863	Harikul Science	HSU012C	1 Mar 23	29 Feb 24
7	SCT Meter	CONDUCTIVITY (umhos/cm)	Horiba	LAQUA-EC210 / HC0J0020	Technology Promotion Association (Thailand~Japan)	23CH1571	14 Dec 23	12 Dec 24
8	Heating Block	CHEMICAL OXYGEN DEMAND	Hanna Instruments Italia Srl.	HI 839800-02 / H 018500 I	Hanna Instruments (Thailand) Ltd.	HIT-2312-0342	10 Mar 23	9 Mar 24
9	Hot Air Oven	TOTAL DISSOLVED SOLIDS TOTAL SUSPENDED SOLIDS	Memmert	UF55 / B212,0411	Technology Promotion Association (Thailand~Japan)	23TM373	11 Apr 23	10 Apr 24
10	Incubator	TOTAL COLIFORM BACTERIA	Binder	KB400 / 20220000022479	DKSH Technology	C31231678	7 Aug 23	6 Aug 24

List of Instrument Certificates for Environmental Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*
11	Inductively Coupled Plasma- Optical Emission Spectrometer(ICP-OES)	BARIUM SODIUM SODIUM ADSORPTION RATIO TITANIUM	Agilent Technologies, USA	5110 VDV(G8015AA) / MY803001	Agilent Technologies (Thailand) Co.,Ltd.	Preventive Maintenance Checklist	13 Nov 23	12 Nov 24
12	pH Meter	pH	Horiba	LAQUA-PH210 / HA1M0043	technology promotion association (thailand-japan)	23CH423	30 Mar 23	28 Mar 24
13	SCT Meter	SALINITY	YSI Environmental	Pro 30 / 18C103132	Technology Promotion Association (Thailand~Japan)	23CH929	26 Jul 23	25 Jul 24
14	Spectrophotometer	CHROMIUM HEXAVALENT	Agilent	Cary 60 G6860A / MY1541009	DQE Services Co.,Ltd.	SP23-021	20 May 23	19 May 24
15	UV-VIS Spectrophotometer	CHEMICAL OXYGEN DEMAND	Hitachi	U-1900 / 2021-064	DQE Services Co.,Ltd.	SP24-008	16 Jan 24	15 Jan 25
16	UV-VIS Spectrophotometer	NITRATE NITROGEN SULPHATE	Hitachi	U-2900 / 21E22-009	DQE Services Co.,Ltd.	SP24-001	4 Jan 24	3 Jan 25
17	Turbidity Meter (Portable)	TURBIDITY (NTU)	Oakton Instruments(China)	T100IR / 1120501017	Technology Promotion Association (Thailand~Japan)	23CH1148	15 Sep 23	13 Sep 24
18	Water Bath	FECAL COLIFORM BACTERIA	Memmert	WNE 14 / L416.0612	Technology Promotion Association (Thailand~Japan)	23TM194	15 Feb 23	14 Feb 24

Due Date of Calibration* : Based on the annual calibration plan. At least 1 time per year.

List of Instrument Certificates for Environmental Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*
1	Atomic Absorption Spectrometer	CADMIUM COPPER IRON LEAD MANGANESE NICKEL ZINC	Agilent Technologies	AA240FS / MY13160001	Agilent Technologies (Thailand) Co.,Ltd.	Preventive Maintenance Checklist	24 Jan 24	23 Jan 25
2	Atomic Absorption Spectrometer	ARSENIC SELENIUM	Perkin Elmer	PinAAde 900F / PFB20031902	Perkin Elmer Co.,Ltd.	WO-02273773	26 Jun 23	25 Jun 24
3	Analytical Balance	FAT OIL AND GREASE	Mettler Toledo	AB204-S/FACT / 1129361010	National Food Institute,Ministry of Industry, Thailand	2303074-001-01	27 May 23	25 May 24
4	Analytical Balance	TOTAL DISSOLVED SOLIDS	Mettler Toledo	XSR205DU / C210685394	Technology Promotion Association (Thailand-Japan)	23MM113	26 Apr 23	25 Apr 24
5	Analytical Balance	TOTAL SUSPENDED SOLIDS	Mettler Toledo	XSR205DU / C009071872	Technology Promotion Association (Thailand-Japan)	23MM112	26 Apr 23	25 Apr 24
6	DO Meter	BIOCHEMICAL OXYGEN DEMAND	YSI	5100 / 11B 101863	Technology Promotion Association (Thailand-Japan)	24TW39	21 Feb 24	20 Feb 25
7	Mercury Analyzer	MERCURY	NIC. Japan	RA-4500 / 17780278	Coax Group Corporation Ltd.	Preventive Maintenance Report	11 Jul 23	10 Jul 24
8	Inductively Coupled Plasma- Optical Emission Spectrometer(ICP-OES)	BARJUM SODIUM ADSORPTION RATIO TITANIUM	Agilent Technologies, USA	5110 VDV(G8015AA) / MY803001	Agilent Technologies (Thailand) Co.,Ltd.	Preventive Maintenance Checklist	13 Nov 23	12 Nov 24
9	pH Meter	pH	EcoSense	pH100A / JC04744	Technology promotion association (thailand-japan)	23CH525	27 Apr 23	26 Apr 24
10	SCT Meter	SALINITY	YSI Environmental	Pro 30 / 23A104805	Technology Promotion Association (Thailand-Japan)	23CH535	27 Apr 23	27 Apr 24
11	SCT Meter	CONDUCTIVITY (umhos/cm)	YSI Environmental	Pro 30 / 18C103132	Technology Promotion Association (Thailand-Japan)	23CH929	26 Jul 23	25 Jul 24

List of Instrument Certificates for Environmental Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*
12	Spectrophotometer	CHROMIUM HEXAVALENT	Agilent	Cary 60 G6860A / MY1541009	DQE Services Co.,Ltd.	SP23-021	20 May 23	19 May 24
13	UV-VIS Spectrophotometer	CHEMICAL OXYGEN DEMAND	Hitachi	U-1900 / 2021-064	DQE Services Co.,Ltd.	SP24-008	16 Jan 24	15 Jan 25
14	UV-VIS Spectrophotometer	NITRATE NITROGEN SULPHATE	Hitachi	U-2900 / 21E22-009	DQE Services Co.,Ltd.	SP24-001	4 Jan 24	3 Jan 25
15	Turbidity Meter (Portable)	TURBIDITY (NTU)	Oakton Instruments(China)	T100IR / 1120501017	Technology Promotion Association (Thailand-Japan)	23CH1148	15 Sep 23	13 Sep 24

Due Date of Calibration* : Based on the annual calibration plan. At least 1 time per year.

List of Instrument Certificates for Environmental Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*
1	Atomic Absorption Spectrometer	CADMIUM CHROMIUM HEXAVALENT CHROMIUM TRIVALENT COPPER IRON LEAD MANGANESE NICKEL ZINC	Agilent Technologies	AA240FS / MY13160001	Agilent Technologies (Thailand) Co.,Ltd.	Preventive Maintenance Checklist	24 Jan 24	23 Jan 25
2	Atomic Absorption Spectrometer	ARSENIC SELENIUM	Perkin Elmer	PinAAcle 900F / PFB20031902	Perkin Elmer Co.,Ltd.	WO-02273773	26 Jun 23	25 Jun 24
3	Analytical Balance	FAT OIL AND GREASE	Mettler Toledo	AB204-S/FACT / 1129361010	National Food Institute,Ministry of Industry, Thailand	2303074-001-01	27 May 23	25 May 24
4	Analytical Balance	TOTAL DISSOLVED SOLIDS	Mettler Toledo	XSR205DU / C210685394	Technology Promotion Association (Thailand~Japan)	23MM113	26 Apr 23	25 Apr 24
5	Analytical Balance	TOTAL SUSPENDED SOLIDS	Mettler Toledo	XSR205DU / C009071872	Technology Promotion Association (Thailand~Japan)	23MM112	26 Apr 23	25 Apr 24
6	Continuous Flow Analyzer(CFA)	CYANIDE AS HCN	Skalar Analytical B.V., the Netherlands	San++5000-02 / 182688	DKSH (Thailand) Ltd.	Service Report/Test Report WO-00018067	20 Feb 24	19 Feb 25
7	DO Meter	BIOCHEMICAL OXYGEN DEMAND	YSI	5100 / 11B 101863	Technology Promotion Association (Thailand~Japan)	24TW39	21 Feb 24	20 Feb 25
8	Digestion Units	TOTAL KJELDAHL NITROGEN	Foss Tecator	2520 Auto / 91794469	National Food Institute Ministry of Industry, Thailand	2302413-001-01	30 Mar 23	29 Mar 24
9	Mercury Analyzer	MERCURY	NIC. Japan	RA-4500 / 17780278	Coax Group Corporation Ltd.	Preventive Maintenance Report	11 Jul 23	10 Jul 24
10	Incubator	TOTAL COLIFORM BACTERIA	Binder	KB400 / 20220000000391	DKSH Technology	C31231210	9 Jun 23	8 Jun 24

List of Instrument Certificates for Environmental Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*
11	Inductively Coupled Plasma- Optical Emission Spectrometer(ICP-OES)	BARIUM SODIUM ADSORPTION RATIO TITANIUM	Agilent Technologies, USA	5110 VDV(G8015AA) / MY803001	Agilent Technologies (Thailand) Co.,Ltd.	Preventive Maintenance Checklist	13 Nov 23	12 Nov 24
12	Kjelttec System Distilling Unit	TOTAL KJELDAHL NITROGEN	Foss Tecator (Labtec)	KT200 / 91790524	FOSS South East Asia	9810	8 Feb 24	7 Feb 25
13	pH Meter	pH	EcoSense	pH100A / JC04744	Technology promotion association (thailand~japan)	23CH525	27 Apr 23	26 Apr 24
14	SCT Meter	SALINITY	YSI Environmental	Pro 30 / 23A104805	Technology Promotion Association (Thailand~Japan)	23CH535	27 Apr 23	27 Apr 24
15	SCT Meter	CONDUCTIVITY (umhos/cm)	YSI Environmental	Pro 30 / 18C103132	Technology Promotion Association (Thailand~Japan)	23CH929	26 Jul 23	25 Jul 24
16	Spectrophotometer	HYDROGEN SULPHIDE	Agilent	Cary 60 G6860A / MY1541009	DQE Services Co.,Ltd.	SP23-021	20 May 23	19 May 24
17	UV-VIS Spectrophotometer	CHEMICAL OXYGEN DEMAND FORMALDEHYDE	Hitachi	U-1900 / 2021-064	DQE Services Co.,Ltd.	SP24-008	16 Jan 24	15 Jan 25
18	UV-VIS Spectrophotometer	NITRATE NITROGEN SULPHATE	Hitachi	U-2900 / 21E22-009	DQE Services Co.,Ltd.	SP24-001	4 Jan 24	3 Jan 25
19	Turbidity Meter (Portable)	TURBIDITY (NTU)	Oakton Instruments(China)	T100IR / 1120501017	Technology Promotion Association (Thailand~Japan)	23CH1148	15 Sep 23	13 Sep 24

Due Date of Calibration* : Based on the annual calibration plan. At least 1 time per year.



Cert.No.: 23CH423
Page.: 1 of 3

Certificate of Calibration

Equipment : pH Meter
Manufacturer : Horiba
Model : LAQUA-PH210
Serial No. : HA1M0043
ID No. : UAE.EFM.013/2565(EFM.pH.03/65)
Condition As-Received: Used Item
Received Date : 28 March 2023
Calibration Date : 30 March 2023
Reference : 2303-1001WSC-6
Submitted by : United Analyst and Engineering Consultant Co.,Ltd
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong, Bangkok 10260
Ambient Temperature : (25 ± 2.5) °C
Relative Humidity : (50 ± 15) %
Calibration Procedure : In - house method :
- CP-CH5 by direct measurement with standard
voltage calibrator and direct measurement with
certified reference material (CRM)
- CP-CH8 by comparison with standard thermometer

Calibrated by : Warakorn Lernagatrakul

Approved by :
Approved Signatory

(/) Malee Butkruea
() Sathip Meangmai
() Warakorn Lernagatrakul

Issue Date : 31 March 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Condition of this calibration result

1. Reference Standard Instrument : -

Instrument	Serial No.	ID No.	Cert. No.	Due Date
1) Document Process Calibrator	54030049	130RC116	22E2769	24 Aug 2023
2) Ref. Standard Thermometer	4982054	110RC044	2211306	27 Oct 2023

This certification is traceable to the International System of Unit maintained at:-
- Traceable to National Institute of Metrology (Thailand), NIMT

2. Certified Reference Materials : The measurement results are traceable to SI through CPA chem Ltd.,
ANSI-ASQ National Accreditation Board, Accredited No. AR-1835

Buffer Solution	Manufacturer	Lot No.	Exp. date
pH 4.008	CPA chem	863832	28 Dec 2024
pH 6.987	CPA chem	826589	09 July 2023
pH 10.010	CPA chem	863835	28 Dec 2023

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration Results

Function : mV Measurement

Performing standard curve by Fluke at pH (4.7)(7,10)

Unit Under Calibration	Nominal Value	Standard Voltage Input	Actual Reading		Uncertainty of Measurement (±mV)	Coverage factor k
	pH	mV	mV	pH		
	pH Meter S/N.: HA1M0043	4.00	177.48	177.4	4.01	0.058
7.00		0.00	-0.2	7.00	0.058	2.00
7.00		0.00	-0.2	7.00	0.058	2.00
10.00		-177.48	-177.6	10.01	0.058	2.00

Malee

เอกสารไม่ควบคุม



Cert.No.: 23CH423
Page.: 3 of 3

Calibration Results

Function : pH Measurement

Performing three buffers standard curve by using buffer nominal pH (4.7)(7,10)

Unit Under Calibration	Standard pH Buffer Solution	Actual pH Reading	Actual mV Reading (mV)	Uncertainty of pH measurement (±)	Coverage factor k
pH Electrode S/N : Q92M0159	4.008	4.01	180.8	0.0085	2.05
	6.987	6.99	6.8	0.011	2.00
	6.987	7.00	6.3	0.011	2.00
	10.010	10.00	-168.9	0.0092	2.00

Function : Temperature Measurement

(*) Without adjustment

This equipment was connected with Temperature Probe;

- Model : 9652-10D
- Serial No. : Q92M0159

Dimension of probe;

- Length : 107 mm
- Diameter : 16 mm
- Immersion Depth : 100 mm

Calibration Point (°C)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty of measurement (± °C)	Coverage factor k
25.0	25.004	25.0	-0.004	0.13	2.00
30.0	30.001	30.0	-0.001	0.13	2.00
35.0	35.004	35.0	-0.004	0.13	2.00

Remark : - UUC* = Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-o0o-

Malee

เอกสารไม่ควบคุม



Cert.No.: 23CH525
Page.: 1 of 3

Certificate of Calibration

Equipment : pH Meter
Manufacturer : EcoSense
Model : pH100A
Serial No. : JC04744
ID No. : UAE.EFM.058/2566(EFM.pH.01/66)
Condition As-Received: Used Item
Received Date : 26 April 2023
Calibration Date : 27 April 2023
Reference : 2304-0707WSC-3
Submitted by : United Analyst and Engineering Consultant Co.,Ltd
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong, Bangkok 10260
Ambient Temperature : (25 ± 2.5) °C
Relative Humidity : (50 ± 15) %
Calibration Procedure : In - house method :
- CP-CH5 by direct measurement with standard
voltage calibrator and direct measurement with
certified reference material (CRM)
- CP-CH8 by comparison with standard thermometer

Calibrated by : Warakorn Lernagatrakul

Approved by :
Approved Signatory

(/) Malee Butkruea
() Sathip Meangmai
() Warakorn Lernagatrakul

Issue Date : 9 May 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Cert.No.: 23CH525
Page.: 2 of 3

Condition of this calibration result

1. Reference Standard Instrument :-

Instrument	Serial No.	ID No.	Cert. No.	Due Date
1) Document Process Calibrator	54030049	130RC116	22E2769	24 Aug 2023
2) Ref. Standard Thermometer	4982054	110RC044	22I1306	27 Oct 2023

This certification is traceable to the International System of Unit maintained at:-
- Traceable to National Institute of Metrology (Thailand), NIMT

2. Certified Reference Materials

The measurement results are traceable to SI through CPA chem Ltd.,
ANSI-ASQ National Accreditation Board, Accredited No. AR-1835

Buffer Solution

Manufacturer

Lot No.

Exp. date

pH 4.008	CPA chem	863832	28 Dec 2024
pH 6.987	CPA chem	826589	09 July 2023
pH 10.010	CPA chem	863835	28 Dec 2023

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration Results

Function : mV Measurement

Performing standard curve by Fluke at pH (4.7)(7,10)

Unit Under Calibration	Nominal Value	Standard Voltage Input	Actual Reading		Uncertainty of Measurement	Coverage factor
	pH	mV	mV	pH	(±mV)	k
pH Meter S/N.: JC04744	4.00	177.48	177	4.01	0.58	2.00
	7.00	0.00	0	7.00	0.58	2.00
	7.00	0.00	0	7.00	0.58	2.00
	10.00	-177.48	-177	10.01	0.58	2.00



Cert.No.: 23CH525
Page.: 3 of 3

Calibration Results

Function : pH Measurement

Performing three buffers standard curve by using buffer nominal pH (4.7)(7,10)

Unit Under Calibration	Standard pH Buffer Solution	Actual pH Reading	Actual mV Reading (mV)	Uncertainty of pH measurement (±)	Coverage factor k
pH Electrode S/N.:230308SIA605377	4.008	4.01	172	0.0079	2.00
	6.987	7.00	-1	0.011	2.00
	6.987	6.99	-1	0.011	2.00
	10.010	10.01	-177	0.0095	2.00

Function : Temperature Measurement

(*) Without adjustment

This equipment was connected with Temperature Probe;

- Model :	-
- Serial No. :	230308SIA605377
- Dimension of probe;	
- Length :	110 mm
- Diameter :	12 mm
- Immersion Depth :	100 mm

Calibration Point (°C)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty of measurement (± °C)	Coverage factor k
25.0	25.002	25.1	0.098	0.13	2.00
30.0	30.003	30.1	0.097	0.13	2.00
35.0	35.002	35.1	0.098	0.13	2.00

Remark : - UUC* = Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k , providing a level of confidence of approximately 95 %.

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL: 0-2717-3000-29 FAX: 0-2719-9484



Cert.No.: 23CH535
Page.: 1 of 3

Certificate of Calibration

Equipment : Conductivity Meter
Manufacturer : YSI
Model : Pro 30
Serial No. : 23A104805
ID No. : UAE EFM.067/2566 (EFM.SCT.03/66)
Condition As-Received: New Item
Received Date : 26 April 2023
Calibration Date : 27 April 2023
Reference : 2304-0706WSC-3
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road, Bangchak,
Phrakhanong, Bangkok 10260
Ambient Temperature : (25 ± 2.5) °C
Relative Humidity : (50 ± 15) %
Calibration Procedure: In-house method :
- CP-CH6 by direct measurement
with certified reference material (CRM)
- CP-CH8 by comparison with standard thermometer

Calibrated by : Walalak Sinthean

Approved by :

Saithip
Approved Signatory

() Malee Butkruea
(✓) Saithip Meangmai
() Warakorn Lerngagtrakul

Issue Date : 3 May 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3: Equipment Calibration and Testing Services.



Cert.No.: 23CH535
Page.: 2 of 3

Condition of this result of calibration

1. Reference Standard Instrument :-

Instrument	Serial No.	ID No.	Certificate No.	Due date
1) Thermometer	1963878	130RC095	22I1140	12 Sep 2023
2) Ref. Std. Thermometer	4982054	110RC044	22I1306	27 Oct 2023

This certification is traceable to the International System of Unit maintained at:-

- Traceable to National Institute of Metrology (Thailand), NIMT

2. Certified Reference Materials :-

- Conductivity calibration solution, CPA chem Ltd., The measurement results are traceable to SI
through CPA chem Ltd., ANSI-ASQ National Accreditation Board, Accredited No. AR-1835

Conductivity Solution	Manufacturer	Lot No.	Exp. date
1413.0 µS/cm	CPA Chem	826595	09 July 2023
12.880 mS/cm	CPA Chem	826596	08 July 2023

- Control Conductivity calibration solution temperature by Water bath (25±0.1) °C
3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration results

Function : Conductivity Measurement

(*) After Adjustment at 1413.0 µS/cm

Conductivity Electrode Serial No.: 23A100616

Standard Conductivity Solution	Before Adjustment UUC* Reading	After Adjustment UUC* Reading	Uncertainty of Measurement (±)	Coverage factor k
1413.0 µS/cm	1443 µS/cm	1413 µS/cm	9.2 µS/cm	2.00
12.880 mS/cm	13.04 mS/cm	12.80 mS/cm	0.086 mS/cm	2.00

Remark : - UUC* = Unit Under Calibration

เอกสารไม่ควบคุม

เอกสารไม่ควบคุม



Cert.No.: 23CH535
Page.: 3 of 3

Calibration Results

Function : Temperature Measurement

(*) Without adjustment

This equipment was connected with Temperature Probe;

- Model : PRO 30 COND-T
- Serial No. : 23A100616

Dimension of probe;

- Length : 94 mm
- Diameter : 2.5 mm
- Immersion Depth : 90 mm

Calibration Point (°C)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty of Measurement (± °C)	Coverage factor k
25.0	25.001	25.0	-0.001	0.13	2.00
30.0	30.001	30.0	-0.001	0.13	2.00
35.0	35.002	35.0	-0.002	0.13	2.00

Remark : - UUC* = Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-000-

Saitip

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/6 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000-29 FAX. 0-2719-9484



Cert.No.: 23CH929
Page.: 1 of 3

Certificate of Calibration

Equipment : Conductivity Meter
Manufacturer : YSI
Model : Pro 30
Serial No. : 18C103132
ID No. : UAE.EFM.204/2561(ENV.SCT.02/61)
Condition As-Received: Used Item
Received Date : 25 July 2023
Calibration Date : 26 July 2023
Reference : 2307-0791WSC-3
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road, Bangchak,
Phrakhanong, Bangkok 10260
Ambient Temperature : (25 ± 2.5) °C
Relative Humidity : (50 ± 15) %
Calibration Procedure: In -house method :
- CP-CH6 by direct measurement
with certified reference material (CRM)
- CP-CH8 by comparison with standard thermometer

Calibrated by : Warakorn Lerngagtrakul

Approved by : 
Approved Signatory

() Malee Bulkruea
(✓) Saitip Meangmai
() Warakorn Lerngagtrakul

Issue Date : 27 July 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Cert.No.: 23CH929
Page.: 2 of 3

Condition of this result of calibration

1. Reference Standard Instrument :-

Instrument	Serial No.	ID No.	Certificate No.	Due date
1) Thermometer	1963878	130RC095	2211140	12 Sep 2023
2) Ref. Std. Thermometer	A5A339	60RC020	2211251	11 Oct 2023

This certification is traceable to the International System of Unit maintained through:-
- Technology Promotion Association (Thailand-Japan)

2. Certified Reference Materials :-

- Conductivity calibration solution, CPA chem Ltd., The measurement results are traceable to SI
through CPA chem Ltd., ANSI-ASQ National Accreditation Board, Accredited No. AR-1835

Conductivity Solution	Manufacturer	Lot No.	Exp. date
1413.0 µS/cm	CPA Chem	885122	28 Mar 2024
12.880 mS/cm	CPA Chem	885123	28 Mar 2024

- Control Conductivity calibration solution temperature by Water bath (25.0 ± 0.1) °C

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration results

Function : Conductivity Measurement

(*) After Adjustment at 1413.0 µS/cm

Conductivity Electrode Serial No.: 18C100559

Standard Conductivity Solution	Before Adjustment UUC* Reading	After Adjustment UUC* Reading	Uncertainty of Measurement (±)	Coverage factor k
1413.0 µS/cm	1176 µS/cm	1414 µS/cm	9.2 µS/cm	2.00
12.880 mS/cm	10.11 mS/cm	12.30 mS/cm	0.086 mS/cm	2.00

Remark : - UUC* = Unit Under Calibration

เอกสารไม่ควบคุม



Cert.No.: 23CH929
Page.: 3 of 3

Calibration Results

Function : Temperature Measurement

(*) Without adjustment

This equipment was connected with Temperature Probe;

- Model : PRO30 COND-T
- Serial No. : 18C100559

Dimension of probe;

- Length : 7 mm
- Diameter : 2.5 mm
- Immersion Depth : 90 mm

Calibration Point (°C)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty of Measurement (± °C)	Coverage factor k
25.0	25.003	25.0	-0.003	0.13	2.00
30.0	30.002	30.0	-0.002	0.13	2.00
35.0	35.002	35.0	-0.002	0.13	2.00

Remark : - UUC* = Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-000-

เอกสารไม่ควบคุม



Cert.No.: 23CH1148
Page.: 1 of 2

Certificate of Calibration

Equipment : Turbidity Meter
Manufacturer : Oakton
Model : T100IR
Serial No. : 1120501017
ID. No. : UAE.WAT.056/2563
Condition As-Received: Used Item
Received Date : 13 September 2023
Calibration Date : 14 September 2023
Reference : 2309-0458DSC-1
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong, Bangkok 10260
Ambient Temperature : (25 ± 2.5) °C
Relative Humidity : (50 ± 20) %
Calibration Procedure : In - house method : CP-CH11
based on direct measurement by
using Formazin standard solution
Calibrated by : Walalak Sirinthean
Approved by :
() Saitthip Meangmai
(✓) Warakorn Lernagatrakul
() Ponpan Paipim
Issue Date : 15 September 2023

The Uncertainties are for a confidence probability of approximately 95%.

This certificate may not be reproduced other than in full, except with the prior written approval of the head of Calibration and Testing Equipment Services.

เอกสารไม่ควบคุม

A 0011853



Cert.No.: 23CH1148
Page.: 2 of 2

Condition of this calibration result

1. Reference Standard Instruments :

This certification is traceable to the International System of unit (SI unit) through:-
- Technology Promotion Association (Thailand-Japan).

Instruments	Serial No.	ID No.	Certificate No.	Due date
1) Thermo-Hygrograph	1103328	130EC010	23C1361	13 June 2024
2) Electronic Balance	1124013382	140RC006	23MM18	20 Feb 2024

2. Standard Material : The Formazin suspension has been prepared gravimetric from

Material	Manufacturer	Lot No.	Assay
1) Hexamethylenetetramine	HIMEDIA	0000493947	99.65%
2) Hydrazinium Sulfate	HIMEDIA	0000522014	99.40%

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration result

Performing five - Formazin suspension standard curve by using 0,20,100,400,800 NTU
Turbidity Meter Serial Number : 1120501017

Standard Formazine suspension (NTU)	UUC* Reading (NTU)	Uncertainty of Measurement (± NTU)	Coverage Factor k
0	0.00	0.0067	2.00
20	20.3	0.39	2.00
100	101	0.76	2.00
400	401	1.5	2.05
800	800	2.1	2.23

Remark - UUC* = Unit Under Calibration
- NTU = Nephelometric Turbidity Units

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-000-

เอกสารไม่ควบคุม

a 1179917



Cert.No.: 23CH1571
Page.: 1 of 3

Certificate of Calibration

Equipment : Conductivity Meter
Manufacturer : Horiba
Model : LAQUA-EC210
Serial No. : HCOJ0020
ID No. : UAE.EFM.078/2564(EFM.SCT.04/64)
Condition As-Received: Used Item
Received Date : 13 December 2023
Calibration Date : 14 December 2023
Reference : 2312-0277WSC-2
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong, Bangkok 10260
Ambient Temperature : (25 ± 2.5) °C
Relative Humidity : (50 ± 15) %
Calibration Procedure : In - house method :
- CP-CH6 by direct measurement
with certified reference material (CRM)
- CP-CH8 by comparison with standard thermometer
Calibrated by : Walalak Sirinthean
Approved by :
(✓) Saitthip Meangmai
() Warakorn Lernagatrakul
() Ponpan Paipim
Issue Date : 18 December 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

A 0061667



Cert.No.: 23CH1571
Page.: 2 of 3

Condition of this result of calibration

1. Reference Standard Instrument :-

Instrument	Serial No.	ID No.	Certificate No.	Due date
1) Thermometer	9549224	130RC003	23I435	10 Apr 2024
2) Ref. Std.Thermometer	4982054	110RC044	23I908	26 July 2024

- This Certification is traceable to SI Through Technology Promotion Association (Thailand - Japan)

2. Certified Reference Materials :-

- Conductivity calibration solution, CPA chem Ltd., The measurement results are traceable to SI through CPA chem Ltd., ANSI-ASQ National Accreditation Board, Accredited No. AR-1835

Conductivity Solution	Manufacturer	Lot No.	Exp. date
1413.0 µS/cm	CPA Chem	931955	30 Sep 2024
12.880 mS/cm	CPA Chem	913597	14 July 2024

- Control Conductivity calibration solution temperature by Water bath (25±0.1) °C

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration results

Function : Conductivity Measurement

(*) After Adjustment at 1413.0 µS/cm

Conductivity Electrode Serial No.: 9B0K0167

Standard Conductivity Solution	Before Adjustment UUC* Reading	After Adjustment UUC* Reading	Uncertainty of Measurement (±)	Coverage factor k
1413.0 µS/cm	1402 µS/cm	1413 µS/cm	9.2 µS/cm	2.00
12.880 mS/cm	12.30 mS/cm	12.64 mS/cm	0.086 mS/cm	2.00

Remark - UUC* = Unit Under Calibration

Saitthip

a 1193974



Cert.No.: 23CH1571
Page.: 3 of 3

Calibration Results

Function : Temperature Measurement

(*) Without adjustment

This equipment was connected with Temperature Probe;

- Model : 9383
- Serial No. : 990K0167

Dimension of probe;

- Length : 104 mm
- Diameter : 16 mm
- Immersion Depth : 90 mm

Calibration Point (°C)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty of Measurement (± °C)	Coverage factor k
25.0	25.001	25.0	-0.001	0.13	2.00
30.0	30.001	30.0	-0.001	0.13	2.00
35.0	35.000	35.1	0.100	0.13	2.00

Remark : - UUC* = Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k , providing a level of confidence of approximately 95 %.

-000-

Saitip

a 1193975



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000-29 FAX. 0-2719-9484



Cert.No.: 23MM112
Page.: 1 of 3

Certificate of Calibration

Equipment : Electronic Balance
Manufacturer : Mettler Toledo
Model : XSR205
Serial No. : C009071872
ID No. : UAE.WAO.012/2563
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phakhanong,
Bangkok 10260
Location : Balance Room
Received order : 26 April 2023
Calibration Date : 26 April 2023
Ambient Temperature : 15 °C to 40 °C
Relative Humidity : 30 % to 90 %
Calibrated by : Man Pattanapongpaiboon
Approved by :
() Ponthippa Tameyakul
() Malee Butkruea
(x) Suwit Imjai
Issue Date : 2 May 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3: Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : Electronic Balance
Condition As-Received : Used Item
Reference : 2304-0459OC-1

Cert.No.: 23MM112
Page: 2 of 3

Procedure used :-

Calibration were conducted using in-house calibration procedure CP-OB01 according to direct measurement method against standard weight.

Condition of this result of calibration

1. Reference standard instruments:-

Instruments	Model	Serial No.	ID No.	Test report No.	Due date
1) Standard Weight Set (E2)	15884	24053	70RC007	MM-0010-22	20 Jan 2024

- This certificate is valid only to the item calibrated on date and place of calibration.
- This result of calibration was made on requested at the point specified by customer.
- This certificate is not certified for any commercial transaction.
- This certification is traceable to the International System of Unit.

Result of calibration () Without Adjustment (x) After Adjustment by Internal Calibration

Range capacity : 0 g to 81 g Resolution 0.00001 g
81 g to 220 g Resolution 0.0001 g

Before Adjustment :

Applied Weight (g)	Balance Reading (g)	Correction (g)	Measurement Uncertainty (± mg)	Coverage Factor (k)
80	80.00005	-0.00005	0.15	2.00
200	199.9999	+0.0001	0.29	2.00

After Adjustment :

1. Determination of the standard deviation of weighing machine (n = 10)

Applied Weight (g)	Standard Deviation of Reading (g)
80	0.000007
200	0.00000

เอกสารไม่ควบคุม




Equipment : Electronic Balance
Condition As-Received : Used Item
Reference : 2304-0459OC-1

Cert.No.: 23MM112
Page: 3 of 3

Result of calibration

2. Effect of off center loading

A mass of 100 g was placed to various position on the pan.
The weighing machine reading error obtained is given in the table

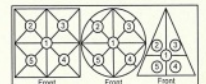

Maximum difference between off-center and central loading (g)

Position 1 (g)	Position 2 (g)	Position 3 (g)	Position 4 (g)	Position 5 (g)
-0.0001	-0.0001	0.0000	-0.0001	-0.0001

3. Departure from nominal value

Applied Weight (g)	Balance Reading (g)	Correction (g)	Measurement Uncertainty (± mg)	Coverage Factor (k)
Unload	0.00000	0.00000	0.014	2.13
0.05	0.05001	-0.00001	0.015	2.09
0.1	0.10001	-0.00001	0.015	2.09
1	1.00001	-0.00001	0.018	2.04
5	5.00003	-0.00003	0.026	2.00
20	20.00006	-0.00006	0.045	2.00
50	50.00006	-0.00006	0.080	2.00
80	80.00004	-0.00004	0.15	2.00
100	100.00000	0.00000	0.16	2.00
150	150.00000	0.00000	0.29	2.00
200	200.00000	0.00000	0.29	2.00

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k , providing a level of confidence of approximately 95 %.

-000-

เอกสารไม่ควบคุม



Cert.No.: 23MM113
Page: 1 of 3

Certificate of Calibration

Equipment : Electronic Balance
Manufacturer : Mettler Toledo
Model : XSR205
Serial No. : C210685394
ID No. : UAE.WAO.0102565
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phakhanong,
Bangkok 10260
Location : Balance Room
Received order : 26 April 2023
Calibration Date : 26 April 2023
Ambient Temperature : 15 °C to 40 °C
Relative Humidity : 30 % to 90 %
Calibrated by : Man Pattanapongpaiboon
Approved by :
() Pornthippa Tameyakul
() Malee Butkruea
(✓) Suwit Imjai
Issue Date : 2 May 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : Electronic Balance
Condition As-Received : Used Item
Reference : 2304-0459OC-2
Procedure used :-

Cert.No.: 23MM113
Page: 2 of 3

Calibration were conducted using in-house calibration procedure CP-OB01 according to direct measurement method against standard weight.

Condition of this result of calibration

1. Reference standard instruments:-

Instruments	Model	Serial No.	ID No.	Test report No.	Due date
1) Standard Weight Set (E2)	15884	24053	70RC007	MM-0010-22	20 Jan 2024

- This certificate is valid only to the item calibrated on date and place of calibration.
- This result of calibration was made on requested at the point specified by customer.
- This certificate is not certified for any commercial transaction.
- This certification is traceable to the International System of Unit.

Result of calibration () Without Adjustment (*) After Adjustment by Internal Calibration

Range capacity :	0 g to 81 g	Resolution	0.00001 g
81 g to 220 g	Resolution	0.0001 g	

Before Adjustment :

Applied Weight	Balance Reading	Correction	Measurement Uncertainty	Coverage Factor
(g)	(g)	(g)	(± mg)	(k)
80	79.99992	+0.00008	0.15	2.00
200	199.9995	+0.0005	0.29	2.00

After Adjustment :

1. Determination of the standard deviation of weighing machine (n = 10)

Applied Weight	Standard Deviation of Reading (g)
(g)	
80	0.000007
200	0.00004

เอกสารไม่ควบคุม




Equipment : Electronic Balance
Condition As-Received : Used Item
Reference : 2304-0459OC-2

Cert.No.: 23MM113
Page: 3 of 3

Result of calibration

2. Effect of off center loading

A mass of 100 g was placed to various position on the pan.
The weighing machine reading error obtained is given in the table

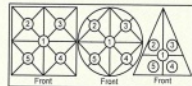

Maximum difference between
off-center and central loading
(g)
0.0001

Position 1 (g)	Position 2 (g)	Position 3 (g)	Position 4 (g)	Position 5 (g)
-0.0001	-0.0001	0.0000	-0.0001	-0.0001

3. Departure from nominal value

Applied Weight	Balance Reading	Correction	Measurement Uncertainty	Coverage Factor
(g)	(g)	(g)	(± mg)	(k)
Unload	0.00000	0.00000	0.014	2.11
0.05	0.04999	+0.00001	0.015	2.09
0.1	0.09999	+0.00001	0.015	2.07
1	1.00000	0.00000	0.018	2.04
5	5.00000	0.00000	0.026	2.00
20	20.00002	-0.00002	0.045	2.00
50	50.00002	-0.00002	0.080	2.00
80	80.00002	-0.00002	0.15	2.00
100	100.00000	0.00000	0.17	2.00
150	150.00000	0.00000	0.29	2.00
200	199.99999	+0.00001	0.29	2.00

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k , providing a level of confidence of approximately 95 %.

-000-

เอกสารไม่ควบคุม



Cert. No.: 23TM194
Page: 1 of 3

Certificate of Calibration

Equipment : Water Bath
Manufacturer : Memmert
Model : WNE 14
Serial No. : L416.0612
ID No. : UAE.MIC.003/2560
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260
Location : Microbiology Laboratory
Received Order : 15 February 2023
Calibration Date : 15 February 2023
Ambient Temperature : (26 ± 10) °C
Relative Humidity : (50 ± 30) %
Calibrated by : Suwit Imjai
Approved by :
() Pornthippa Tameyakul
(✓) Malee Butkruea
Issue Date : 24 February 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : Water Bath
Condition As-Received : Used Item
Reference : 2302-0295OC-3
Procedure Used :-

Cert. No.: 23TM194
Page : 2 of 3

Calibration were conducted using in-house calibration procedure CP-OT04 according to direct measurement method with Data Acquisition which connected with Industrial Platinum Resistance Thermometer (IPRT).

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument Model Serial No. Cert. No. Due Date
1) Data Acquisition 34972A MY59003411 22LM165 26 Nov 2023


2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

	Environmental		AC Voltage Supply
	(°C)	(%R.H.)	(Volt)
Beginning of Calibration	22	65	231
Finished of Calibration	22	63	230

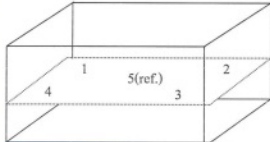

Front

Position :	Ref. Std. ID No.:
1	4804539-001
2	4804539-002
3	4804539-003
4	4804539-004
5(ref.)	4804539-005

เอกสารไม่ควบคุม



Equipment : Water Bath
Condition As-Received : Used Item
Reference : 2302-0295OC-3
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source

Cert. No.: 23TM194
Page : 3 of 3

Calibration point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Average* Standard Reading (°C)				
			Position				
			1	2	3	4	5 (ref.)
44.5	44.5	44.6	44.520	44.509	44.498	44.552	44.530

Calibration point (°C)	Uniformity (°C)	Stability (± °C)	Uncertainty (± °C)	Coverage Factor k
44.5	0.077	0.037	0.15	2

Average* : The average of 30 values in each position.

Uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Stability : One-half of the greatest maximum difference of measured temperature at any one probe.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity.

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k , providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000-29 FAX. 0-2719-9484



Cert. No.: 23TM373
Page : 1 of 3

Certificate of Calibration

Equipment : Hot Air Oven
Manufacturer : Memmert
Model : UF 55
Serial No. : B212.0411
ID No. : UAE.WAO.005/2556
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260
Location : Lab Floor 2
Received Order : 11 April 2023
Calibration Date : 11 - 12 April 2023
Ambient Temperature : (26 ± 10) °C
Relative Humidity : (50 ± 30) %
Calibrated by : Krisda Malee

Approved by :
() Pornthippa Tameyakul
() Malee Butkruea
() Suwit Imjai

Issue Date : 24 April 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3: Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม

A 0053359



Equipment : Hot Air Oven
Condition As-Received : Used Item
Reference : 2304-0156OC-1
Procedure Used :-

Cert. No.: 23TM373
Page : 2 of 3

Calibration were conducted using calibration procedure CP-OT02 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD) and Thermocouple Type T.

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument Model Serial No. Cert. No. Due Date
1) Data Acquisition 34972A MY59003411 22LM165 26 Nov 2023


2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Close

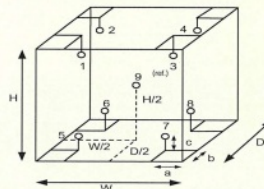

Probe Installation Details : Dimension of Chamber :
a = 5.0 cm D = 0.50 m
b = 5.0 cm W = 0.80 m
c = 5.0 cm H = 0.75 m
Capacity = 0.30 m³

Environment during calibration		
	Beginning	Finished
Temp. (°C)	27	28
REL.Humid. (%)	45	44
AC Supply (Volt)	221	220

Ref. Std. ID No.: @ Calibration Point		
Position :	(120 to 180) °C	(104) °C
1	18-20TC-01	20RTD-2/1
2	18-20TC-02	20RTD-2/2
3	18-20TC-03	20RTD-2/3
4	18-20TC-04	20RTD-2/4
5	18-20TC-05	20RTD-2/5
6	18-20TC-06	20RTD-2/6
7	18-20TC-07	20RTD-2/7
8	18-20TC-08	20RTD-2/8
9 (ref.)	18-20TC-09	20RTD-2/9

เอกสารไม่ควบคุม

a 1158261



Equipment : Hot Air Oven
Condition As-Received : Used Item
Reference : 2304-0156OC-1
Result of Calibration : (°) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Close

Cert. No.: 23TM373
Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor k
104.0	104.0	104.0	0.054	0.59	0.95	2
120.0	120.0	120.0	0.12	0.89	1.5	2
180.0	180.0	180.0	0.12	1.5	2.5	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	1	2	3	4	5	6	7	8	9 (ref.)	
104.0	104.512	104.016	104.542	104.407	103.704	103.729	104.167	104.158	104.001	0.42
120.0	120.317	119.768	120.524	120.232	119.363	119.209	119.888	119.797	119.735	1.1
180.0	180.878	179.819	181.357	180.871	179.303	179.139	180.230	180.055	179.960	1.1

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-00-

เอกสารไม่ควบคุม
a 1158260



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3 : EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000 FAX. 0-2719-9484

Cert.No.: 23TW90
Page.: 1 of 2

Certificate of Testing

Equipment : DO Meter
Manufacturer : Horiba
Model : LAQUA-DO210
Serial No. : HE1M0020
ID No. : UAE.EFM.019/2565(EFM.DO.03/65)
Received Date : 28 March 2023
Test Date : 29 March 2023
Reference : 2303-0996WSC-7
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road, Bangchak,
Phrakhanong, Bangkok 10260
Laboratory Condition : Temperature (25 ± 5) °C
Humidity (50 ± 20) %
Test Procedure : In - house method : CP-CH9
by Comparison Technique with Azide Modification Method
Tested by : Walalak Sirinthean

Approved by :
Approved Signatory

(/) Malee Butkrua
() Saitip Meangmai
() Warakorn Lernagatrakul

Issue Date : 31 March 2023

เอกสารไม่ควบคุม
a 1158260



Cert.No.: 23TW90
Page.: 2 of 2

Condition of this result of calibration

- Reference Standard Instruments :
This certification is traceable to the International System of Unit through the reference standards laboratory of Industrial Calibration Center, Technology Promotion Association (Thailand-Japan).

Instruments	Serial No.	ID No.	Certificate No.	Due Date
1) Burette	-	130BU10	23CG1172	22 Mar 2025
2) Balance	1126143764	140RC004	22MM50	20 Sep 2023

2. Standard Material :-

Material	Manufacturer	Lot.No.	Assay
Sodium Thiosulfate pentahydrate	Merck	AM1763316	100.2%

Result : Dissolved Oxygen Meter Adjustment With Air 100 %

Dissolved Oxygen Probe No.: 9K1K0080

Titration Method (Azide Modification Method) (mg/L)	DO Meter Reading (mg/L)	Standard Deviation (mg/L)
8.16	8.15	0.0055

This report was certified only for the instrument we tested. It is allowable to use for study the system efficiency. The environmental impact control and present to organization it may concerned. Intend to use for advertising and referral purpose is prohibited. This report may not be reproduced other in full, without written approval of the laboratory.

-00-

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3 : EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000 FAX. 0-2719-9484



Cert. No.: 23LM81
Page.: 1 of 2

Certificate of Calibration

Equipment : DO Meter With Sensor
Manufacturer : Horiba
Model : LAQUA-DO210
Serial No. : HE1M0020
ID No. : UAE.EFM.019/2565(EFM.DO.03/65)
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260
Location : TPA On Site Calibration Laboratory
Received Order : 28 March 2023
Calibrated Date : 31 March 2023
Ambient Temperature : (26 ± 10) °C
Relative Humidity : (50 ± 30) %
AC Line Voltage : (220 ± 22) V

Calibrated by : Malee Butkrua

Approved by :
Approved Signatory

() Pornthippa Tameyakul
(/) Suwit Imjai

Issue Date : 12 April 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม
a 1158260



Equipment : DO Meter With Sensor
Condition As-Received : Used Item
Reference : 2303-0996WSC-8
Procedure Used :-

Cert. No.: 23LM81
Page.: 2 of 2

Calibration were conducted using in-house calibration procedure CP-OT01 according to comparison with Industrial Platinum Resistance Thermometer (IPRT) into Temperature Bath.
The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Model	Serial No.	Cert. No.	Due Date
1) Digital Thermometer	1502A	A52847	2211325	31 Oct 2023

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

Result of Calibration :- (°) Without Adjustment

Function : Temperature measurement.

This instrument was connected with temperature sensor, S/N.: 9K1K0080

Calibration Point (°C)	Immersion Depth (mm)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty (± °C)	Coverage Factor k
25.0	90	24.999	25.0	0.001	0.16	2.00
30.0	90	29.997	30.0	0.003	0.16	2.00
35.0	90	35.008	35.0	-0.008	0.16	2.00

UUC* : Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

เจต

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3 : EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000 FAX. 0-2719-9484

Cert.No.: 24TW39
Page.: 1 of 2

Certificate of Testing

Equipment : DO Meter
Manufacturer : YSI
Model : 5100
Serial No. : 11B 101863
ID No. : UAE.WAO.004/2554
Received Date : 20 February 2024
Test Date : 21 February 2024
Reference : 2402-06290SC-1
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udonsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260
Laboratory Condition : Temperature (25 ± 5) °C
Humidity (50 ± 20) %
Test Procedure : In - house method : CP-CH9
by Comparison Technique with Azide Modification Method
Tested by : Walalak Sirithuan
Approved by :
Approved Signatory
() Pornthippa Tameyakul
() Unnoppol Harachai
(✓) Sathip Meangmai
Issue Date : 22 February 2024

เอกสารไม่ควบคุม



Cert.No.: 24TW39
Page.: 2 of 2

Condition of this result of calibration

1. Reference Standard Instruments :

This certification is traceable to the International System of Unit through the reference standards laboratory of Industrial Calibration Center, Technology Promotion Association (Thailand-Japan).

Instruments	Serial No.	ID No.	Certificate No.	Due Date
1. Burette	-	130BU10	23CG1172	22 Mar 2025
2. Balance	14233821	110RC001	23MM405	16 July 2024

2. Standard Material :-

Material	Manufacturer	Lot.No.	Assay
Sodium Thiosulfate pentahydrate	Merck	AM1763316	100.2%

Result : Dissolved Oxygen Meter Adjustment With Air 100 %
Dissolved Oxygen Probe No.: 22B100125

Titration Method (Azide Modification Method) (mg/L)	DO Meter Reading (mg/L)	Standard Deviation (mg/L)
8.20	8.19	0.0055

This report was certified only for the instrument we tested. It is allowable to use for study Intend to use for advertising and referral purpose is prohibited. This report may not be reproduced other in full, without written approval of the laboratory

-o0o-

เอกสารไม่ควบคุม

FOSS

Customer Service Report

Date: 9 Feb 2024		Report No: 9810																																									
Customer: UAE		Address: Bangkok																																									
Instrument: KT200		Serial: 91790524																																									
Hours Start	Travel To Customer 08:00	Labour 09:30	Travel From Customer 16:30																																								
Finish	09:30	17:00	18:00																																								
<table><thead><tr><th colspan="2">Application</th><th colspan="2">Special</th><th colspan="2">Job Type</th><th colspan="2">Standard</th></tr></thead><tbody><tr><td>Normal</td><td>x</td><td>Courtesy Visit</td><td>x</td><td>Installation</td><td>x</td><td>Training</td><td>x</td></tr><tr><td>Distributor</td><td>x</td><td>PMA Onboarding</td><td>x</td><td>Quote</td><td>x</td><td>In House</td><td>x</td></tr><tr><td>Internal</td><td>x</td><td>Warranty</td><td>x</td><td>Repair</td><td>x</td><td>PM</td><td>x</td></tr><tr><td>Digital Service</td><td>x</td><td>Sales Support</td><td>x</td><td>Remote</td><td>x</td><td>Other</td><td>x</td></tr></tbody></table>				Application		Special		Job Type		Standard		Normal	x	Courtesy Visit	x	Installation	x	Training	x	Distributor	x	PMA Onboarding	x	Quote	x	In House	x	Internal	x	Warranty	x	Repair	x	PM	x	Digital Service	x	Sales Support	x	Remote	x	Other	x
Application		Special		Job Type		Standard																																					
Normal	x	Courtesy Visit	x	Installation	x	Training	x																																				
Distributor	x	PMA Onboarding	x	Quote	x	In House	x																																				
Internal	x	Warranty	x	Repair	x	PM	x																																				
Digital Service	x	Sales Support	x	Remote	x	Other	x																																				
PO/Quote Number: If applicable																																											
PMA Type: FOSSCAT		Contract No. If applicable																																									
<table><thead><tr><th>Details of Work / Test</th><th>Condition / Status</th></tr></thead><tbody><tr><td># PM KT200</td><td></td></tr><tr><td>- ตรวจสอบและปรับตั้ง DO meter</td><td></td></tr><tr><td>- ตรวจสอบและปรับตั้ง DO meter 100 mL</td><td></td></tr><tr><td>- Calibrate 30 mL - 20 mL</td><td></td></tr><tr><td>- ตรวจสอบ PM kit</td><td></td></tr><tr><td>- ตรวจสอบ PM kit</td><td></td></tr><tr><td># ตรวจสอบ SERIAL Model ตรวจสอบและปรับตั้ง DO meter</td><td></td></tr><tr><td>10000725 SERIAL Model complete 1 PC</td><td></td></tr></tbody></table>				Details of Work / Test	Condition / Status	# PM KT200		- ตรวจสอบและปรับตั้ง DO meter		- ตรวจสอบและปรับตั้ง DO meter 100 mL		- Calibrate 30 mL - 20 mL		- ตรวจสอบ PM kit		- ตรวจสอบ PM kit		# ตรวจสอบ SERIAL Model ตรวจสอบและปรับตั้ง DO meter		10000725 SERIAL Model complete 1 PC																							
Details of Work / Test	Condition / Status																																										
# PM KT200																																											
- ตรวจสอบและปรับตั้ง DO meter																																											
- ตรวจสอบและปรับตั้ง DO meter 100 mL																																											
- Calibrate 30 mL - 20 mL																																											
- ตรวจสอบ PM kit																																											
- ตรวจสอบ PM kit																																											
# ตรวจสอบ SERIAL Model ตรวจสอบและปรับตั้ง DO meter																																											
10000725 SERIAL Model complete 1 PC																																											
Instrument Ready for Use <input checked="" type="checkbox"/> OK <input type="checkbox"/> Not OK If not OK - Comment																																											
Part No:	Batch	Description	Qty																																								
10009865	14.12.2020	Foss PM kit KT200 beclot Analyser 100	1																																								
I confirm this report is accurate and complete																																											
Signed FOSS	Signed Customer																																										
Name	Name																																										
Would you be willing to participate in a brief survey in order to tell us how we performed?																																											
Email																																											

เอกสารไม่ควบคุม

Calibration Certificate

Certificate No.: 2303074-001-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road,
Bangchack, Prakanong, Bangkok 10260

Page 1 of 3

Equipment: Electronic Balance
Manufacturer: METTLER TOLEDO
Model: AB204-S/FACT
Serial No.: 1129361010
ID No.: UAE.WAS.002/2552
Order No.: 2303074
Operation No.: 2303074-001
Date of Receipt: 26 May 2023
Date of Calibration: 26 May 2023

Calibrated by Mr.Pheraphat Tuanjit
Scientist
Approved by *P. Jengcharit*
(Miss Preeyaporn Jaengkarnkit)
Vice President, Department of Laboratory Services
Responsible for the Technical Management Team
Date of Issue: 29 May 2023

The uncertainties are for a confidence probability of approximately 95%

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

FCS-009 Revision: 01 Date: 20-04-65

Calibration Report

Certificate No.: 2303074-001-01
Equipment: Electronic Balance
Model: AB204-S/FACT
Serial No.: 1129361010
Capacity: 220 g
Manufacturer: METTLER TOLEDO
Resolution: 0.0001 g
ID No.: UAE.WAS.002/2552

Page 2 of 3

Date of Calibration: 26 May 2023
Environment Condition: Ambient Temperature: 23.7 ± 0.1 °C Relative Humidity: 61 ± 2.2 %
Place of Calibration: Room 108 Balance Room, UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.

Condition of Equipment: Good Condition

Condition of This Results of Calibration:

1. Calibration Method: NFI Method W-MA-001 In-House Method based on UKAS Lab 14 : 2019

2. Reference Standards:

Reference Standard	Model	Serial No.	Calibrated By	Certificate No.	Due Date
Standard Weight Class E2	1mg to 200g	8505567572	TCS	M23040535	8 April 2024
Instrument	Model	Serial No.	Calibrated By	Certificate No.	Due Date
Thermo-Hygro Meter	608-H1	NFI.BTH 018/23	Quality Reborn	QR23-0491	21 February 2024

3. This certification is traceable to SI UNIT

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

Calibration Results:


1. Repeatability of Reading:

Nominal Value (g)	Standard Deviation of Reading (g)
100	0.000048
200	0.000048

2. Off-Center Error:

A mass of 100 g was placed and moved to various position on pan.

The balance reading obtained is given in the table.

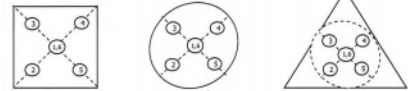

1	2	3	4	5	6	(Maximum Difference)
(g)	(g)	(g)	(g)	(g)	(g)	(g)
99.9996	99.9995	99.9995	99.9999	99.9999	99.9997	0.0003

FCS-012 Revision: 01 Date: 20-04-65

Calibration Report

Certificate No.: 2303074-001-01
Equipment: Electronic Balance
Model: AB204-S/FACT
Serial No.: 1129361010
Capacity: 220 g
Manufacturer: METTLER TOLEDO
Resolution: 0.0001 g
ID No.: UAE.WAS.002/2552

Date of Calibration: 26 May 2023

Page 3 of 3

Calibration Results: (Continued)

Calibration Range: 0-200 g

Calibration Adjustment: Internal Calibration

3. Departure from Nominal Value:

Nominal Value (g)	Standard Value (g)	Average Reading (g)	Correction (g)	Uncertainty (±g)	Coverage Factor k
Unload	0.00000	0.0000	0.0000	0.000088	2.00
0.01	0.01000	0.0100	0.0000	0.000088	2.00
0.05	0.05000	0.0500	0.0000	0.000088	2.00
0.1	0.10001	0.0999	0.0001	0.000088	2.00
0.2	0.20001	0.1999	0.0001	0.000088	2.00
0.5	0.50002	0.5000	0.0000	0.000088	2.00
1	1.00000	1.0000	0.0000	0.000089	2.00
2	2.00002	2.0000	0.0000	0.000089	2.00
5	5.00002	5.0000	0.0000	0.000090	2.00
10	10.00001	9.9999	0.0001	0.000091	2.00
20	20.00003	20.0000	0.0000	0.000095	2.00
50	50.00003	49.9999	0.0001	0.00011	2.00
70	70.00006	69.9999	0.0002	0.00013	2.00
100	100.00006	99.9999	0.0002	0.00016	2.00
150	150.00009	149.9999	0.0002	0.00021	2.00
200	200.00016	199.9998	0.0004	0.00028	2.00

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

----- End -----

FCS-012 Revision: 01 Date: 20-04-65

Certificate of Calibration

Equipment: Incubator
Model: KB 400
Serial No.(or ID): 20220000000301
Manufacturer: Binder
Condition: In Condition
Shelves(pc.): 5
Certificate No.: C31231210
Issued Date: 10 June 2023
Job No.: KSPR2308771
Page: 1 of 3
Ventilation Valve: None

Customer: United Analyst and Engineering Consultant Company Limited.
3 Soi Udomsuk 41 Sukhumvit Road,
Bangkak, Prakanong, Bangkok 10260 Thailand.

Environment Condition: Temperature: 21 °C ± 0.7 °C
Humidity: 58 %RH ± 4.0 %RH
Voltage: 229 VAC ± 1.3 VAC

Calibration Place: United Analyst and Engineering Consultant Company Limited.
3 Soi Udomsuk 41 Sukhumvit Road, (Microbiology Laboratory)
Bangkak, Prakanong, Bangkok 10260 Thailand.

Calibration By: Mr. Amornthep Phumppo

Calibration Date: 09 June 2023

The Method used: In house method, CAL-WI-16, base on TLAS-G20

Traceability: This certificate is traceable to the SI Units maintained by National Institute of Metrology (NIMT), Thailand through DKSH Technology Limited.
Certificate No. C10220023

(Mr. Amornthep Phumppo)

Person in charge

(Mr. Udon Srichana)

Authorized signatory

This certificate is issued the units of measurement according to the International System of Units (SI). It provides traceability of measurement to international or national standard or other recognized national standard laboratories.
The measurement uncertainty stated is the expanded uncertainty which is obtained from the standard uncertainty multiplied by the coverage factor (k=2) to provide a level of confidence of approximately 95%. It is determined in accordance with the Guide to Expression of Uncertainty in Measurement (GUM).


These results may be affected by deviations from specified conditions. The results relate only to the items tested, calibrated or sampled. The report shall not be reproduced except in full without approval of DKSH Technology Limited.

บริษัท ดีเคเอส อีเซีย จำกัด
DKSH Technology Limited
2533 ซอยสุขุมวิท 41 แขวงคลองเตย เขตคลองเตย กรุงเทพมหานคร 10260
Phone: +66 (0) 2422 8688 Fax: +66 (0) 2422 8545
Email: info.certificate@dksh.com Website: www.dksh.com/scientific-thailand

Delivering Growth - in Asia and Beyond.

เอกสารไม่ควบคุม

CAL-FM-C31-10: 12 Sep 2022

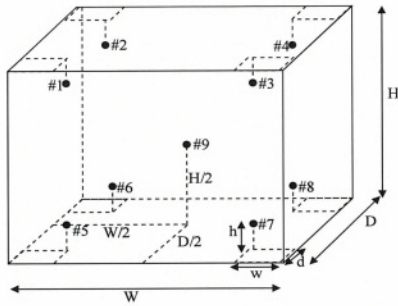

Standard Installation Locations

Volume (Calibration Zone)= 188 (Liters)

Inside chamber: W = 65 (cm) D = 48 (cm) H = 127 (cm)

Standard Locations (#1, #2, #3, #4): w = 7 (cm) d = 5 (cm) h = 15 (cm)

Standard Locations (#5, #6, #7, #8): w = 7 (cm) d = 5 (cm) h = 15 (cm)

#9: Geometric center of the chamber

Position of Std	#1	#2	#3	#4	#5	#6	#7	#8	#9
Channel of Logger	1	2	3	4	5	6	7	8	9

Definitions

Indicating Temperature: The average reading of indicating device which forms the integral part of the enclosure.**Measured Temperature:** The average reading of standards at any positions or location.**Measured Uniformity:** The maximum difference of measured temperatures between of any probes and the measured temperature at the reference location which are observed at same time or at close observation time as possible to determine the temperature pattern or homogeneity with the chamber at steady-state. The reference probe is preferably located in the geometric center of the chamber.**Measured Stability:** The one-half of greatest maximum difference of measured temperatures at any one probe.**Overall Variation:** The difference of maximum and minimum measured temperatures throughout observation time.

บริษัท ดีเคเอส อีเซีย จำกัด
DKSH Technology Limited
2533 สุขุมวิท ถนน, กรุงเทพฯ, ประเทศไทย 10260
2533 Sukhumvit Road, Bangkok, Prathanon, Bangkok 10260
Phone: +66 2639 7000 Email: info.calibration@dksh.com Website: www.dksh.com/scientific-thailand

Delivering Growth - in Asia and Beyond.

เอกสารไม่ควบคุม

CAL-FM-C31-10: 12 Sep 2022

Calibration Results:

Without adjustment

Measurement Temperature at Spread Locations, Indicating of Unit Under Calibration: 35.0 °C

Locations	Measured Temperature (°C)	Correction of UUC (°C)	Uncertainty (± °C)
#1	35.23	0.23	0.23
#2	35.14	0.14	0.23
#3	35.13	0.13	0.23
#4	34.95	-0.05	0.23
#5	35.08	0.08	0.23
#6	35.05	0.05	0.23
#7	34.96	-0.04	0.23
#8	34.79	-0.21	0.23
#9	35.00	0.00	0.23

Temperature Distribution

Desired (°C)	Setting (°C)	Indicating (°C)	Measured Temperature at Spread Locations (°C)									Uncertainty (± °C)
35.0	35.0	35.0	#1	#2	#3	#4	#5	#6	#7	#8	#9	
			35.23	35.14	35.13	34.95	35.08	35.05	34.96	34.79	35.00	0.23

Chamber Characterization

Indicating (°C)	Measured Uniformity (°C)	Measured Stability (± °C)	Overall Variation (°C)
35.0	0.26	0.06	0.48

Note: * Maximum uncertainty of the each position

The End of Certificate

บริษัท ดีเคเอส อีเซีย จำกัด
DKSH Technology Limited
2533 สุขุมวิท ถนน, กรุงเทพฯ, ประเทศไทย 10260
2533 Sukhumvit Road, Bangkok, Prathanon, Bangkok 10260
Phone: +66 2639 7000 Email: info.calibration@dksh.com Website: www.dksh.com/scientific-thailand

Delivering Growth - in Asia and Beyond.

เอกสารไม่ควบคุม

CAL-FM-C31-10: 12 Sep 2022



Certificate of Calibration

Equipment: Cooled Incubator
Model: KB 400
Serial No.(or ID): 2022000022479
Manufacturer: Binder
Condition: New
Shelves(pc.): 5

Certificate No.: C31231678
Issued Date: 10 August 2023
Job No.: WO-00002652
Page: 1 of 3
Ventilation Valve: None

Customer: United Analyst and Engineering Consultant Company Limited.
3 Soi Udomsuk 41 Sukhumvit Road,
Bangkok, Prakanong, Bangkok 10260 Thailand.

Environment Condition: Temperature: 25 °C ± 1.9 °C
Humidity: 49 %RH ± 5.3 %RH
Voltage: 232 VAC ± 1.2 VAC

Calibration Place: United Analyst and Engineering Consultant Company Limited. (Control Area)
3 Soi Udomsuk 41 Sukhumvit Road,
Bangkok, Prakanong, Bangkok 10260 Thailand.

Calibration By: Mr. Thanakrit Raksapol
Calibration Date: 07 August 2023
The Method used: In house method, CAL-WI-16, base on TLAS-G20
Traceability: This certificate is traceable to the SI Units maintained by National Institute of Metrology (NIMT), Thailand through DKSH Technology Limited.
Certificate No. C10230019

Thakrit
(Mr. Thanakrit Raksapol)
Person in charge

Udon
(Mr. Udon Srichana)
Authorized signatory


This certificate is issued the units of measurement according to the International System of Units (SI). It provides traceability of measurement to international or national standard or other recognized national standard laboratories.
The measurement uncertainty stated is the expanded uncertainty which is obtained from the standard uncertainty multiplied by the coverage factor (k=2) to provide a level of confidence of approximately 95%. It is determined in accordance with the Guide to Expression of Uncertainty in Measurement (GUM).
These results may be affected by deviations from specified conditions. The results relate only to the items tested, calibrated or sampled. The report shall not be reproduced except in full without approval of DKSH Technology Limited.

บริษัท ดีเคเอส อีเซีย จำกัด
DKSH Technology Limited
2533 สุขุมวิท ถนน, กรุงเทพฯ, ประเทศไทย 10260
2533 Sukhumvit Road, Bangkok, Prathanon, Bangkok 10260
Phone: +66 2639 7000 Email: info.calibration@dksh.com Website: www.dksh.com/scientific-thailand

Delivering Growth - in Asia and Beyond.

เอกสารไม่ควบคุม

CAL-FM-C31-10: 12 Sep 2022

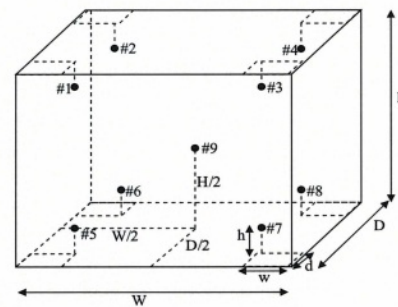

Standard Installation Locations

Volume (Calibration Zone)= 193 (Liters)

Inside chamber: W = 65 (cm) D = 49 (cm) H = 127 (cm)

Standard Locations (#1, #2, #3, #4): w = 7 (cm) d = 5 (cm) h = 15 (cm)

Standard Locations (#5, #6, #7, #8): w = 7 (cm) d = 5 (cm) h = 15 (cm)

#9: Geometric center of the chamber

Position of Std	#1	#2	#3	#4	#5	#6	#7	#8	#9
Channel of Logger	301	302	303	304	305	306	307	308	309

Definitions

Indicating Temperature: The average reading of indicating device which forms the integral part of the enclosure.**Measured Temperature:** The average reading of standards at any positions or location.**Measured Uniformity:** The maximum difference of measured temperatures between of any probes and the measured temperature at the reference location which are observed at same time or at close observation time as possible to determine the temperature pattern or homogeneity with the chamber at steady-state. The reference probe is preferably located in the geometric center of the chamber.**Measured Stability:** The one-half of greatest maximum difference of measured temperatures at any one probe.**Overall Variation:** The difference of maximum and minimum measured temperatures throughout observation time.

บริษัท ดีเคเอส อีเซีย จำกัด
DKSH Technology Limited
2533 สุขุมวิท ถนน, กรุงเทพฯ, ประเทศไทย 10260
2533 Sukhumvit Road, Bangkok, Prathanon, Bangkok 10260
Phone: +66 2639 7000 Email: info.calibration@dksh.com Website: www.dksh.com/scientific-thailand

Delivering Growth - in Asia and Beyond.

เอกสารไม่ควบคุม

CAL-FM-C31-10: 12 Sep 2022

Calibration Results:

Without adjustment

Measurement Temperature at Spread Locations, Indicating of Unit Under Calibration: 35.0 °C

Locations	Measured Temperature (°C)	Correction of UUC. (°C)	Uncertainty (± °C)
#1	35.11	0.11	0.23
#2	35.04	0.04	0.23
#3	35.03	0.03	0.23
#4	35.13	0.13	0.23
#5	35.02	0.02	0.23
#6	35.07	0.07	0.23
#7	34.97	-0.03	0.23
#8	34.97	-0.03	0.23
#9	35.10	0.10	0.23

Temperature Distribution

Desired (°C)	Setting (°C)	Indicating (°C)	Measured Temperature at Spread Locations (°C)									Uncertainty (± °C)*
			#1	#2	#3	#4	#5	#6	#7	#8	#9	
35.0	35.0	35.0	35.11	35.04	35.03	35.13	35.02	35.07	34.97	34.97	35.10	0.23

Chamber Characterization

Indicating (°C)	Measured Uniformity (°C)	Measured Stability (± °C)	Overall Variation (°C)
35.0	0.16	0.04	0.22

Note: * Maximum uncertainty of the each position

The End of Certificate

Statements of conformity:

This conformity certificate documents the validity of the following statements of conformity based on the measurement results of corresponding calibration certificate:

The correction of indication determined during calibration are under given measurement and environmental conditions and considering the expanded measurement uncertainty (coverage probability 95%) within the specification. The given measurement uncertainty already includes other all effects by according to the standard method, TLAS-G20. Therefore, those parameters have not been assessed separately.

Tolerance and Decision rules:

Assessment of the conformity of the measurement device are done based on direct comparison of the relevant measurement results with the tolerances and decision rule are prescribed by the customer.

Decision rule : ☐ Choice A Binary Statement for Simple Acceptance Rule (w = 0), Specific Risk < 50% PFA.

☒ Choice B Non-binary statement with guard band (w = 1 U), Pass or Fail Specific Risk < 2.5% PFA and Condition Pass or Condition Fail Specific Risk < 50% PFA.

☐ Choice C Customer defined, Customers may define arbitrary multiple of r to have applied as guard band (w = r U) .
; PFA - Probability of False Accept

(Mr. Udon Srichana)
Authorized signatory

Without adjustment

Desired Temperature : 35.0°C Tolerances : 0.5 °C

Measurement Temperature at Spread Locations, Indicating of Unit Under Calibration: 35.0 °C

Locations	Measured (°C)	Correction* (°C)	Guard band (W) (± °C)	Tolerance (± °C)	Conformity
#1	35.11	0.11	0.23	0.5	Pass
#2	35.04	0.04	0.23	0.5	Pass
#3	35.03	0.03	0.23	0.5	Pass
#4	35.13	0.13	0.23	0.5	Pass
#5	35.02	0.02	0.23	0.5	Pass
#6	35.07	0.07	0.23	0.5	Pass
#7	34.97	-0.03	0.23	0.5	Pass
#8	34.97	-0.03	0.23	0.5	Pass
#9	35.10	0.10	0.23	0.5	Pass

Correction* = Measured Temperature - Desired Temperature

The validity of the statements of conformity cannot be guaranteed for different places of use, environmental conditions or improper use.

The End of Statements of Conformity

บริษัท ดีเคเอส อีเซีย จำกัด
DKSH Technology Limited
2533 สุขุมวิท ถนน แขวงคลองเตย เขตคลองเตย กรุงเทพมหานคร 10260
2533 Sukhumvit Road, Bangkok, Prachinong, Bangkok 10260
Phone: +66 2639 7000 Email: info.calibration@dksh.com Website: www.dksh.com/scientific-thailand

Delivering Growth - In Asia and Beyond.

เอกสารไม่ควบคุม



CAL-FM-C31-10: 12 Sep 2022

บริษัท ดีเคเอส อีเซีย จำกัด
DKSH Technology Limited
2533 สุขุมวิท ถนน แขวงคลองเตย เขตคลองเตย กรุงเทพมหานคร 10260
2533 Sukhumvit Road, Bangkok, Prachinong, Bangkok 10260
Phone: +66 2639 7000 Email: info.calibration@dksh.com Website: www.dksh.com/scientific-thailand

Delivering Growth - In Asia and Beyond.

เอกสารไม่ควบคุม

CAL-FM-C31-10: 12 Sep 2022

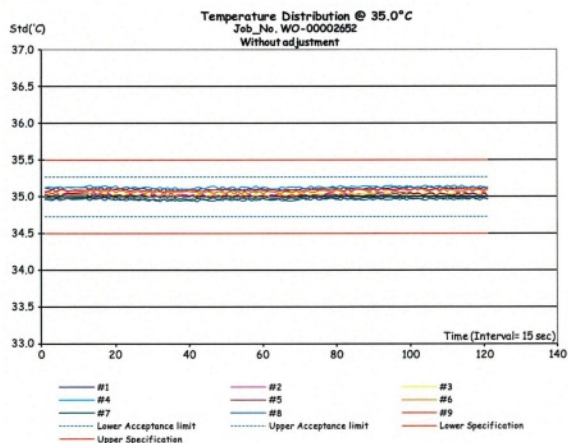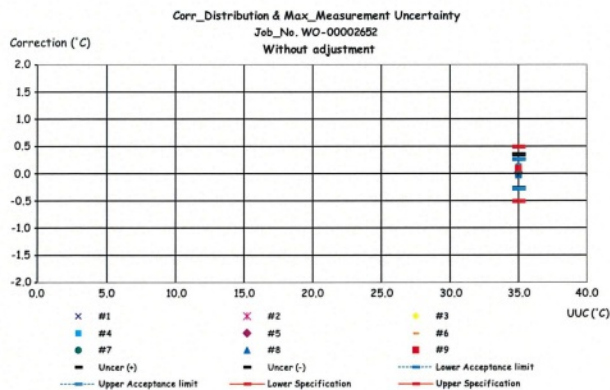

ใบตรวจสอบสภาพเครื่องควบคุมอุณหภูมิ

เลขที่ใบงาน: WO-00002652

ชนิดเครื่องมือ: Cooled Incubator

รุ่น: KB 400

หมายเลขเครื่อง: 20220000022479.000

ตรวจสอบ (รับ)		ตรวจสอบ (ส่ง)		หมายเหตุ
07 Aug 2023		07 Aug 2023		
ปกติ	ไม่ปกติ	ปกติ	ไม่ปกติ	
General				
<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input type="checkbox"/>	1. สายไฟ
<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input type="checkbox"/>	2. การทำงาน Main Switch
<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input type="checkbox"/>	3. การทำงาน Selector Key
<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input type="checkbox"/>	4. การแสดงผล Display
<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input type="checkbox"/>	5. การทำงาน พัดลม
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	6. สภาพ Lever of Ventilation valve
<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input type="checkbox"/>	7. สภาพ Lever door open / close
<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input type="checkbox"/>	8. สภาพ Door seal
<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input type="checkbox"/>	9. การทำงานของระบบ Safety
<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input type="checkbox"/>	10. การทำงานของระบบทำความเย็น
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	11. การทำงานของระบบทำความร้อน
<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input type="checkbox"/>	12. สภาพตู้เครื่อง
<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input type="checkbox"/>	13. สภาพแวดล้อม ณ สถานที่ตั้งเครื่อง

ข้อแนะนำ:

Mr. Thanakrit Raksapol
Service Engineer

เอกสารไม่ควบคุม

บริษัท ดีเคเอส อีเซีย จำกัด
DKSH Technology Limited
2533 สุขุมวิท ถนน แขวงคลองเตย เขตคลองเตย กรุงเทพมหานคร 10260
2533 Sukhumvit Road, Bangkok, Prachinong, Bangkok 10260
Phone: +66 2639 7000 Email: info.calibration@dksh.com Website: www.dksh.com/scientific-thailand

Delivering Growth - In Asia and Beyond.

เอกสารไม่ควบคุม

Certificate No. : HIT-2312-0342

Page : 1 of 2

CERTIFICATE OF CALIBRATION

Equipment : COD Test Tube Heater
Meter Model : HI839800-02 Serial No. : H0185001
Tube Heater : 25 Vial Capacity Accuracy : $\pm 2^{\circ}\text{C}$
Temperature Range : -10°C to 160°C Temperature of Reaction : 150°C
Ambient Temperature : $(25 \pm 2)^{\circ}\text{C}$ Relative Humidity : $(50 \pm 15)\%$ RH
Manufacturer : Hanna Instruments Made in : Romania
Condition As-Received : Used Product Reference : RE230392
Customer name : United Analyst and Engineering Consultant Co., Ltd.
3 Soi Udomsuk 41, Sukhumvit Rd., Bangchak,
Phrakhanong, Bangkok 10260
Received date : 8 March 2023
Calibrate date : 10 March 2023
Issue date : 20 March 2023
Calibrated Location : Hanna Instruments (Thailand) Ltd.
Calibration Procedure : This calibrator was conducted by using in-house: calibration procedure
CP-04 by using certified reference material.

Mr. Pichit Petthong

Calibrated by : ☐ Mr. Jakkapob Pentisan
☐ Mr. Channarong Soinak

Approved by : Mr. Anan Suwanchaisakul

Authorized Signatory
HANNA
instruments
(Thailand) Limited

This certificate was certified only for the instrument we calibrated.

This result of calibration was found accurate on date and place of calibration only.

** This certificate may not be reproduced other than in full, except with the prior written
approval of the head of Hanna Instrument (Thailand).

เอกสารไม่ควบคุม

Condition of this calibration result

Reference Standard Instruments:

Instruments	Model	Serial No.	Certificate No.	Traceable
Data Acquisition Switch Unit	34970A	MY44065265	WK2207-065-1	WK Electric Co., Ltd.

Calibration Result:

Measurement Temperature Source Accuracy for COD Reactor

Capacity (Vial)	Nominal Value ($^{\circ}\text{C}$)	Average Value ($^{\circ}\text{C}$)	\pm Uncertainty ($^{\circ}\text{C}$)	\pm Tolerance of UUC ($^{\circ}\text{C}$)	Acceptance Criteria
25 Vial	150.0	150.3	0.59	2	Pass

Figure: Shows the location of the temperature source.

(1A)	(2A)	(3A)	(4A)	(5A)
149.78 $^{\circ}\text{C}$	150.31 $^{\circ}\text{C}$	150.63 $^{\circ}\text{C}$	149.93 $^{\circ}\text{C}$	150.31 $^{\circ}\text{C}$
(1B)	(2B)	(3B)	(4B)	(5B)
150.35 $^{\circ}\text{C}$	150.18 $^{\circ}\text{C}$	149.93 $^{\circ}\text{C}$	150.18 $^{\circ}\text{C}$	150.21 $^{\circ}\text{C}$
(1C)	(2C)	(3C)	(4C)	(5C)
150.24 $^{\circ}\text{C}$	151.10 $^{\circ}\text{C}$	150.80 $^{\circ}\text{C}$	150.36 $^{\circ}\text{C}$	150.86 $^{\circ}\text{C}$
(1D)	(2D)	(3D)	(4D)	(5D)
150.16 $^{\circ}\text{C}$	149.77 $^{\circ}\text{C}$	150.22 $^{\circ}\text{C}$	150.67 $^{\circ}\text{C}$	150.43 $^{\circ}\text{C}$
(1E)	(2E)	(3E)	(4E)	(5E)
149.94 $^{\circ}\text{C}$	150.44 $^{\circ}\text{C}$	150.06 $^{\circ}\text{C}$	150.63 $^{\circ}\text{C}$	149.29 $^{\circ}\text{C}$

Remark: The Acceptance criteria is the error value plus or minus the Measurement Uncertainty, and then Not
More than the Tolerance value of UUC, therefore concluded that pass.

The reported uncertainty of measurement was based on standard uncertainty multiplied by a coverage factor
 $k=2$, providing a level of confidence of approximately 95%.

** End of certificate **

เอกสารไม่ควบคุม



CERT.No.: HS-U012C

Calibration Date : 1 Mar 23
Submitted by : United Analyst and Engineering Consultant Co., Ltd.
3 Soi Udomsuk 41, Sukhumvit Road, Bangchak,
Phrakhanong, Bangkok.(Head office)

Avg Room Temp : 20°C
Avg Water Temp : 20°C
Air Pressure : 760.00 mmHg
Salinity : 0 ppt

Model : YSI 5100
S/N : 11B101863
Probe : YSI 5010
S/N : 22B100125
ID NO. : -
Air Temp ref : S/N. E00522
Barometric ref : S/N. E00522
Water Temp ref : S/N. 11431
Technician : Kittipong M.

Harikul Science Co.,Ltd.

694 Soi Ratchadapisek 24, Pracharatbamphen,
Samsen-nok, Huaykwang, Bangkok 10310
Tel: 0-2274-2456 Fax: 0-2274-2443
Email: info@harikul.com www.harikul.com

Certificate of Calibration

Calibration Details

Calibration Point	100% air sat. (@20 $^{\circ}\text{C}$, DO = 9.09 mg/l)	(status)	(status)
Measurement 1 (mg/l)	9.09	(PASS)	-
Measurement 2 (mg/l)	9.09	(PASS)	-
Measurement 3 (mg/l)	9.09	(PASS)	-
Measurement 4 (mg/l)	9.09	(PASS)	-
Measurement 5 (mg/l)	9.09	(PASS)	-
Measurement 6 (mg/l)	9.09	(PASS)	-
Measurement 7 (mg/l)	9.08	(PASS)	-
Measurement 8 (mg/l)	9.09	(PASS)	-
Measurement 9 (mg/l)	9.08	(PASS)	-
Measurement 10 (mg/l)	9.09	(PASS)	-

Mean Measurement 9.09 mg/l
Inaccuracy 0.00 mg/l

Overall Status (PASS)

Manufacturer Specification

Accuracy = ± 0.02 mg/l

- This certificate is issued based on the result that are found as shown on date and place of test only.
- The calibration procedure followed in accordance with Harikul Science Co., Ltd.
- This result shall not be used for advertising purpose.

Technician Signature
(Kittipong Maekwong)

Supreecha S.
Laboratory Manager
(Supreecha Sumartam)

เอกสารไม่ควบคุม



Request No. 25-66 / 0323

MTC. ACL.No. 387 / 66

CALIBRATION CERTIFICATE

NOMENCLATURE : 1. Atomic Absorption Spectrophotometer "Agilent Technologies"
Model AA240FS, Serial No. MY13160001
2. Working standard solution "Inorganic Ventures"
Multi Analyte Custom Grade Solution, Lot No. S2-ME2708640
SUBMITTED BY : United Analyst and Engineering Consultant Co., Ltd.
3 Soi Udomsuk41, Sukhumvit Road, Bangchak, Prakanong, Bangkok 10260

CALIBRATION PROCEDURE : 1. Performance Verification of Atomic Absorption Spectrophotometer
(WI-500-02-30)

2. Estimation Uncertainty of Measurement in Analytical Chemistry (QP-513)

CALIBRATION RANGE: 0.02,0.10,0.30,0.50,0.70 mg/l at 228.8 nm.Cd, 0.10,0.20,0.30,0.50,0.70 mg/l at 357.9 nm.Cr,
0.05,0.10,0.30,0.50,0.70 mg/l at 324.7 nm.Cu, 0.10,0.30,0.50,0.70,1.00 mg/l at 248.3 nm.Fe, 0.20,0.50,0.70,1.00,1.50 mg/l
at 217.0 nm.Pb, 0.05,0.10,0.30,0.50,0.70 mg/l at 279.5 nm.Mn, 0.10,0.30,0.50,0.70,1.00 mg/l at 232.0 nm.Ni,
0.05,0.10,0.30,0.50,0.70 mg/l at 213.9 nm.Zn

CALIBRATION DATE : 2 February 2023

REFERENCE MATERIAL : Traceable to NIST "Carlo Erba", "PanReac AppliChem"

Cadmium Lot No. 1152457, Chromium Lot No. 1793249, Copper Batch No. T117098A, Iron Batch No. T126087A,
Lead Lot No. 1227873, Manganese Batch No. T109228A, Nickel Batch No. T270178A, Zinc Batch No. T820140A

AMBIENT CONDITIONS : Temperature 22°C Relative humidity 58 %

The Atomic Absorption Spectrophotometer has been calibrated against Reference
Material traceable to National Institute of Standards and Technology (NIST) by The Analytical Chemistry
Laboratory. The results are attached herewith.

Calibrated by : Danai Srithongkum
(Mr. Danai Srithongkum)

Approved by : Sumittra Deavong
(Miss Sumittra Deavong)
Senior Technical Officer
Acting Director of Analytical Chemistry Laboratory
Ref. 2015266012600366001
Issued Date : 15 February 2023

The results relate only to the items tested/calibrated or value assigned.
Advertising the Report/Certificate and publicity of the results except in full are prohibited unless written permission is obtained from the governor of TISTR.

Head Office

35 Mu 3 Tambon Khlong Ha, Amphoe Khlong Luang,
Changwat Pathumthani 12120, Thailand
Tel. (66) 0 2577 9000
Tel. (66) 0 2577 9009
E-mail : numpai@tistr.or.th Website:www.tistr.or.th

Office/Laboratory

Soi 1C, Bangpoo Industrial Estate, Sukhumvit Road,
Amphoe Muang, Changwat Samutprakan 10280, Thailand
Tel. (66) 0 2323 1672-80 ext. 115, 116
Fax. (66) 0 2323 9165
E-mail : mtc@tistr.or.th

Office

196 Phahonyothin Road, Chatuchak, Bangkok 10900,
Thailand
Tel. (66) 0 2577 8392
Fax. (66) 0 2577 8392
E-mail : sumalee@tistr.or.th

FMBLMTCC02 Rev.4

เอกสารไม่ควบคุม



Request No. 25-66 / 0323

1 / 5

MTC. ACL. No. 387 / 66

CALIBRATION DATA

1. Noise Level

Element	Cd	Cr	Cu	Fe	Pb	Mn	Ni	Zn
Absorbance	0.0020	0.0000	0.0008	0.0000	-0.0009	0.0021	-0.0016	-0.0022
	0.0015	0.0006	0.0005	-0.0009	-0.0014	0.0018	0.0002	-0.0023
	0.0014	0.0006	0.0010	-0.0009	0.0015	0.0008	-0.0004	-0.0015
	0.0021	-0.0008	0.0013	-0.0010	0.0005	0.0005	-0.0008	-0.0004
	0.0020	-0.0012	0.0004	0.0003	-0.0004	0.0001	-0.0024	-0.001
	0.0021	-0.0011	0.0011	0.0003	0.0006	0.0009	-0.0002	-0.0013
	0.0017	-0.0009	0.0001	-0.0015	0.0010	0.0007	0.0001	-0.0016
	0.0024	-0.0012	0.0004	-0.0002	0.0008	-0.0005	-0.0012	-0.0019
	0.0011	-0.0002	0.0015	-0.0004	0.0004	0.0008	-0.0003	-0.0017
	0.0017	0.0000	0.0009	0.0004	0.0001	0.0015	-0.0009	-0.0024
	0.0019	-0.0004	0.0004	0.0000	0.0006	0.0010	-0.0005	-0.0016
	0.0016	-0.0025	0.0003	0.0005	0.0009	-0.0004	-0.0013	-0.0016
	0.0018	-0.0014	0.0001	-0.0009	-0.0006	0.0010	-0.0004	-0.0017
	0.0019	-0.0006	0.0011	-0.0008	0.0011	0.0004	-0.0003	-0.0005
	0.0024	0.0003	0.0005	-0.0012	-0.0002	0.0012	-0.0006	-0.0011
	0.0023	-0.0012	0.0006	-0.0007	0.0002	0.0014	-0.0012	-0.0013
	0.0020	-0.0014	0.0009	-0.0018	0.0003	0.0012	-0.0012	-0.0013
	0.0010	-0.0015	0.0002	0.0004	0.0017	0.0011	-0.0018	-0.0013
	0.0016	-0.0011	0.0013	0.0003	0.0007	0.0026	-0.0006	-0.0006
	0.0001	-0.0007	0.0009	-0.0003	0.0008	0.0008	0.0000	-0.0001
	Average Absorbance	0.002	-0.001	0.001	0.000	0.000	-0.001	-0.001

Continue 2 / 5

INDUSTRIAL METROLOGY AND TESTING SERVICE CENTRE

The results relate only to the items tested/calibrated or value assigned.

Advertising the Report/Certificate and publicity of the results except in full are prohibited unless written permission is obtained from the governor of TISTR.

FMLBL.MTC.002 Rev.4

Head Office
35 Mu 3 Tambon Khlong Ha, Amphoe Khlong Luang,
Changwat Pathumthani 12120, Thailand
Tel. (66) 0 2577 9000
Fax. (66) 0 2577 9009
E-mail : rumpa@tistr.or.th Website:www.tistr.or.th

Office/Laboratory
Soi 1C, Bangpoo Industrial Estate, Sukhumvit Road,
Amphoe Muang, Changwat Samutprakan 10280, Thailand
Tel. (66) 0 2323 1672-80 ext. 115, 116
Fax. (66) 0 2323 9165
E-mail : mtg@tistr.or.th

Office
196 Phahonyothin Road, Chatuchak, Bangkok 10900,
Thailand
Tel. (66) 0 2379 6592
Fax. (66) 0 2379 6592
E-mail : sumalee@tistr.or.th

เอกสารไม่ควบคุม



Request No. 25-66 / 0323

3 / 5

MTC. ACL. No. 387 / 66

3. Trueness

3.1 Reading on wavelength- Cadmium(Cd) at 228.8 nm.

Element	Standard Value of RM (mg/l)	Reading (mg/l)	Error of Measurement (mg/l)	Error of Measurement (%)	Uncertainty (mg/l)
Cd	0.02002	0.021	0.001	4.90	± 0.005
	0.30030	0.298	-0.002	0.77	± 0.005
	0.70070	0.675	-0.026	3.67	± 0.008

3.2 Reading on wavelength- Chromium (Cr) at 357.9 nm.

Element	Standard Value of RM (mg/l)	Reading (mg/l)	Error of Measurement (mg/l)	Error of Measurement (%)	Uncertainty (mg/l)
Cr	0.1001	0.101	0.001	0.90	± 0.009
	0.3003	0.293	-0.007	2.43	± 0.012
	0.7007	0.648	-0.053	7.52	± 0.023

3.3 Reading on wavelength- Copper (Cu) at 324.7 nm.

Element	Standard Value of RM (mg/l)	Reading (mg/l)	Error of Measurement (mg/l)	Error of Measurement (%)	Uncertainty (mg/l)
Cu	0.050	0.046	-0.004	8.00	± 0.003
	0.300	0.289	-0.011	3.67	± 0.009
	0.700	0.674	-0.026	3.71	± 0.020

Continue 4 / 5

INDUSTRIAL METROLOGY AND TESTING SERVICE CENTRE

The results relate only to the items tested/calibrated or value assigned.

Advertising the Report/Certificate and publicity of the results except in full are prohibited unless written permission is obtained from the governor of TISTR.

Head Office
35 Mu 3 Tambon Khlong Ha, Amphoe Khlong Luang,
Changwat Pathumthani 12120, Thailand
Tel. (66) 0 2577 9000
Fax. (66) 0 2577 9009
E-mail : rumpa@tistr.or.th Website:www.tistr.or.th

Office/Laboratory
Soi 1C, Bangpoo Industrial Estate, Sukhumvit Road,
Amphoe Muang, Changwat Samutprakan 10280, Thailand
Tel. (66) 0 2323 1672-80 ext. 115, 116
Fax. (66) 0 2323 9165
E-mail : mtg@tistr.or.th

Office
196 Phahonyothin Road, Chatuchak, Bangkok 10900,
Thailand
Tel. (66) 0 2379 6592
Fax. (66) 0 2379 6592
E-mail : sumalee@tistr.or.th

เอกสารไม่ควบคุม



Request No. 25-66 / 0323

2 / 5

MTC. ACL. No. 387 / 66

2. Precision

Element	Conc. (mg/l)	Absorbance										Ave. Abs.	SD	%RSD
Cd	0.02	0.0085	0.0084	0.0090	0.0089	0.0089	0.0090	0.0086	0.0092	0.0090	0.0089	0.009	0.0003	2.88
	0.30	0.0993	0.1001	0.1007	0.1004	0.1004	0.0995	0.0997	0.0998	0.0999	0.0996	0.100	0.0005	0.45
	0.70	0.2238	0.2229	0.2244	0.2249	0.2243	0.2233	0.2235	0.2231	0.2251	0.2240	0.224	0.0007	0.33
Cr	0.10	0.0088	0.0087	0.0094	0.0086	0.0086	0.0091	0.0099	0.0095	0.0076	0.0085	0.009	0.0006	7.25
	0.30	0.0257	0.0265	0.0255	0.0270	0.0266	0.0258	0.0261	0.0262	0.0274	0.0262	0.026	0.0006	2.25
	0.70	0.0573	0.0590	0.0580	0.0576	0.0578	0.0579	0.0593	0.0599	0.0586	0.0594	0.058	0.0009	1.51
Cu	0.05	0.0083	0.0084	0.0084	0.0075	0.0086	0.0086	0.0081	0.0080	0.0087	0.0092	0.008	0.0005	5.45
	0.30	0.0430	0.0444	0.0426	0.0429	0.0435	0.0432	0.0428	0.0441	0.0427	0.0436	0.043	0.0006	1.41
	0.70	0.0981	0.0992	0.0990	0.0997	0.0977	0.0986	0.0990	0.0982	0.0988	0.0980	0.099	0.0006	0.63
Fe	0.10	0.0109	0.0104	0.0087	0.0100	0.0087	0.0094	0.0102	0.0092	0.0094	0.0100	0.010	0.0007	7.53
	0.50	0.0456	0.0442	0.0450	0.0444	0.0450	0.0455	0.0445	0.0441	0.0446	0.0444	0.045	0.0006	1.27
	1.00	0.0904	0.0901	0.0891	0.0876	0.0873	0.0901	0.0876	0.0886	0.0879	0.0901	0.089	0.0001	1.38
Pb	0.20	0.0093	0.0099	0.0104	0.0102	0.0104	0.0109	0.0102	0.0103	0.0115	0.0117	0.010	0.0007	6.85
	0.70	0.0344	0.0336	0.0336	0.0328	0.0338	0.0346	0.0336	0.0331	0.0343	0.0350	0.034	0.0007	2.02
	1.50	0.0709	0.0718	0.0706	0.0713	0.0698	0.0718	0.0712	0.0713	0.0715	0.0719	0.071	0.0006	0.90
Mn	0.05	0.0115	0.0130	0.0131	0.0127	0.0135	0.0136	0.0124	0.0133	0.0124	0.0130	0.013	0.0006	4.88
	0.30	0.0709	0.0700	0.0714	0.0704	0.0700	0.0705	0.0714	0.0698	0.0694	0.0700	0.070	0.0007	0.96
	0.70	0.1619	0.1633	0.1646	0.1638	0.1646	0.1614	0.1632	0.1614	0.1636	0.1652	0.163	0.0014	0.83
Ni	0.10	0.0113	0.0105	0.0113	0.0114	0.0110	0.0113	0.0117	0.0112	0.0107	0.0117	0.011	0.0004	3.45
	0.50	0.0509	0.0517	0.0508	0.0502	0.0517	0.0516	0.0516	0.0523	0.0518	0.0503	0.051	0.0007	1.36
	1.00	0.0997	0.1006	0.1006	0.1006	0.0996	0.0998	0.1007	0.1000	0.1013	0.0999	0.100	0.0006	0.55
Zn	0.05	0.0315	0.0309	0.0322	0.0304	0.0329	0.0312	0.0313	0.0319	0.0308	0.0311	0.031	0.0007	2.35
	0.30	0.1705	0.1728	0.1688	0.1693	0.1711	0.1704	0.1704	0.1707	0.1708	0.1688	0.170	0.0012	0.70
	0.70	0.3559	0.3572	0.3548	0.3560	0.3559	0.3550	0.3579	0.3552	0.3574	0.3573	0.356	0.0011	0.31

Continue 3 / 5

INDUSTRIAL METROLOGY AND TESTING SERVICE CENTRE

The results relate only to the items tested/calibrated or value assigned.

Advertising the Report/Certificate and publicity of the results except in full are prohibited unless written permission is obtained from the governor of TISTR.

FMLBL.MTC.002 Rev.4

Head Office
35 Mu 3 Tambon Khlong Ha, Amphoe Khlong Luang,
Changwat Pathumthani 12120, Thailand
Tel. (66) 0 2577 9000
Fax. (66) 0 2577 9009
E-mail : rumpa@tistr.or.th Website:www.tistr.or.th

Office/Laboratory
Soi 1C, Bangpoo Industrial Estate, Sukhumvit Road,
Amphoe Muang, Changwat Samutprakan 10280, Thailand
Tel. (66) 0 2323 1672-80 ext. 115, 116
Fax. (66) 0 2323 9165
E-mail : mtg@tistr.or.th

Office
196 Phahonyothin Road, Chatuchak, Bangkok 10900,
Thailand
Tel. (66) 0 2379 6592
Fax. (66) 0 2379 6592
E-mail : sumalee@tistr.or.th

เอกสารไม่ควบคุม



Request No. 25-66 / 0323

4 / 5

MTC. ACL. No. 387 / 66

3.4 Reading on wavelength- Iron (Fe) at 248.3 nm.

Element	Standard Value of RM (mg/l)	Reading (mg/l)	Error of Measurement (mg/l)	Error of Measurement (%)	Uncertainty (mg/l)
Fe	0.100	0.095	-0.005	5.00	± 0.014
	0.500	0.474	-0.026	5.20	± 0.016
	1.000	0.950	-0.050	5.00	± 0.029

3.5 Reading on wavelength- Lead (Pb) at 217.0 nm.

Element	Standard Value of RM (mg/l)	Reading (mg/l)	Error of Measurement (mg/l)	Error of Measurement (%)	Uncertainty (mg/l)
Pb	0.200	0.207	0.007	3.50	± 0.014
	0.700	0.673	-0.027	3.86	± 0.030
	1.500	1.417	-0.083	5.53	± 0.061

3.6 Reading on wavelength- Manganese (Mn) at 279.5 nm.

Element	Standard Value of RM (mg/l)	Reading (mg/l)	Error of Measurement (mg/l)	Error of Measurement (%)	Uncertainty (mg/l)
Mn	0.04995	0.046	-0.004	7.91	± 0.005
	0.29970	0.294	-0.0057	1.90	± 0.007
	0.69930	0.694	-0.0053	0.76	± 0.014

Continue 5 / 5

INDUSTRIAL METROLOGY AND TESTING SERVICE CENTRE

The results relate only to the items tested/calibrated or value assigned.

Advertising the Report/Certificate and publicity of the results except in full are prohibited unless written permission is obtained from the governor of TISTR.

FMLBL.MTC.002 Rev.4

Head Office
35 Mu 3 Tambon Khlong Ha, Amphoe Khlong Luang,
Changwat Pathumthani 12120, Thailand
Tel. (66) 0 2577 9000
Fax. (66) 0 2577 9009
E-mail : rumpa@tistr.or.th Website:www.tistr.or.th

Office/Laboratory
Soi 1C, Bangpoo Industrial Estate, Sukhumvit Road,
Amphoe Muang, Changwat Samutprakan 10280, Thailand
Tel. (66) 0 2323 1672-80 ext. 115, 116
Fax. (66) 0 2323 9165
E-mail : mtg@tistr.or.th

Office
196 Phahonyothin Road, Chatuchak, Bangkok 10900,
Thailand
Tel. (66) 0 2379 6592
Fax. (66) 0 2379 6592
E-mail : sumalee@tistr.or.th

เอกสารไม่ควบคุม



Request No. 25-66 / 0323

5 / 5

MTC. ACL. No. 387 / 66

3.7 Reading on wavelength- Nickel (Ni) at 232.0 nm.

Element	Standard Value of RM (mg/l)	Reading (mg/l)	Error of Measurement (mg/l)	Error of Measurement (%)	Uncertainty (mg/l)
Ni	0.1001	0.103	0.003	2.90	± 0.013
	0.5005	0.501	0.001	0.10	± 0.018
	1.0010	0.987	-0.014	1.40	± 0.032

3.8 Reading on wavelength- Zinc (Zn) at 213.9 nm.

Element	Standard Value of RM (mg/l)	Reading (mg/l)	Error of Measurement (mg/l)	Error of Measurement (%)	Uncertainty (mg/l)
Zn	0.050	0.046	-0.004	8.00	± 0.013
	0.300	0.311	0.011	3.67	± 0.013
	0.700	0.665	-0.035	5.00	± 0.019

Remark : The reported uncertainty is an expanded uncertainty calculated using a coverage factor of 2 (k = 2) which gives a level of confidence of approximately 95%

Calibrated by 1. Dani Srithongkum
(Mr. Danai Srithongkum)
2. Atipat
(Mr. Atipat Ratana)

Approved by Sutadde Deawong
(Miss Sutadde Deawong)
Senior Technical Officer
Acting Director of
Analytical Chemistry Laboratory
Issued Date : 15 February 2023

INDUSTRIAL METROLOGY AND TESTING SERVICE CENTRE
End of Certificate

The results relate only to the items tested/calibrated or value assigned.
Advertising the Report/Certificate and publicity of the results except in full are prohibited unless written permission is obtained from the governor of TISTR.

FM.BL.MTC.002 Rev.4

Head Office
35 Mu 3 Tambon Khlong Ha, Amphoe Khlong Luang,
Changwat Pathumthani 12120, Thailand
Tel. (66) 0 2577 9000
Fax. (66) 0 2577 9009
E-mail : numpai@tistr.or.th Website: www.tistr.or.th

Office/Laboratory
Sri 1C, Bangpoo Industrial Estate, Sukhumvit Road,
Amphoe Muang, Changwat Samutprakan 10280, Thailand
Tel. (66) 0 2323 1672-80 ext. 115, 116
Fax. (66) 0 2323 9165
E-mail : mtrc@tistr.or.th

Office
196 Phahonyothin Road, Chatuchak, Bangkok 10900,
Thailand
Tel. (66) 0 2319 0392
Fax. (66) 0 2319 0392
E-mail : sumalee@tistr.or.th

เอกสารไม่ควบคุม

Agilent
CrossLab
From Insight to Outcome

Agilent 55 240 280 Series Atomic Absorption Spectroscopy Systems

Preventive Maintenance Checklist

Agilent Preventive Maintenance provides factory recommended service for your analytical systems to assure reliable operation and the accuracy of your results.

Delivered by highly trained and certified service engineers using genuine Agilent parts and supplies, Agilent Preventive Maintenance provides everything you need to reduce unplanned downtime and keep your systems operating at their peak. This checklist will be completed at the end of the service and provided to you as a record of the installation.

Note: While non-current production AA instrument and or accessory models are not covered specifically in this document it can be used as a basic reference.

For more information about Agilent Technologies services please visit our web site using the following URL <http://www.agilent.com/en-us/services>


Introduction

Customer Information

- Customers should provide all necessary operating supplies upon request of the engineer.
- A customer representative should be available to the engineer while performing the preventive maintenance procedures.
- Any parts, not included in the Parts Lists section of this document, are not part of the recommended Preventive Maintenance service, nor are they included in the price of this service.
- If a system requires the use of extra or special procedures and/or parts for the maintenance service, then these must be ordered separately and charged as a repair, which may incur additional costs.

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021

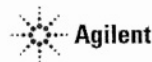

เอกสารไม่ควบคุม

Agilent
CrossLab
From Insight to Outcome

Instrument Preventive Maintenance Checklist

Important Customer Web Links

- For more information about Agilent Technologies services, please visit our website using the following URL: <http://www.agilent.com/en-us/products/crosslab-instrument-services/service-repair>
- To access Agilent University, visit <http://www.agilent.com/crosslab/university/> to learn about training options, which include online, classroom and onsite delivery. A training specialist can work directly with you to help determine your best options.
- A useful Agilent Resource Center web page is available, which includes short videos on maintenance, quick lists of consumables for new instruments, and other valuable information. Check out the Resource Page here: <https://www.agilent.com/en-us/agilentresources>
- Need technical support, FAQs, supplies? – visit our Support Home page at <http://www.agilent.com/search/support>
- Get answers. Share insights. Build connections. Join the Agilent Community at <https://community.agilent.com/welcome>


Service Engineer's Responsibilities

- Contact the customer and ensure that all necessary supplies are available before the preventive maintenance visit.
- Confirm the ability of the instrument to deliver continued safe operation as established via the Agilent AA safe operation flow chart. (Refer directly to the AA 55/240/280 Preventive Maintenance Scope of Work to make this decision.)
- Only select those pages that relate to the system or module being serviced.
- Complete empty fields with the relevant information.
- Complete the relevant checkboxes in the checklist using either a "X" or tick mark "✓".
- Check "Section not applicable" check boxes to indicate services/tasks not delivered, as appropriate.
- Complete the Preventive Maintenance service in the order of the tasks listed.
- Complete the Service Review section together with the customer.
- Complete the fields for page numbers at the foot of each selected page
- Complete the total number of pages field in the Service Completion section
- Ask the customer to sign the Service Completion section including the customer's and your signature.

This information is subject to change without notice.

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021

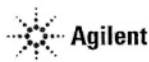

เอกสารไม่ควบคุม

Agilent
CrossLab
From Insight to Outcome

Instrument Preventive Maintenance Checklist

Instrument Maintenance

System Information

☐ Check this box if an instrument configuration report is attached instead of completing the table.

Instrument System Name and ID	
Instrument System Site and Location	UNITED ANALYST AND ENGINEERING CONSULTANT / 2nd Lab FI

List System Component Product Numbers	List the Serial Numbers of each Component
1. G 8432 A	W 0160001
2.	
3.	
4.	
5.	
6.	
7.	
8.	
9.	

Preparation, Safe operation and Initial performance checks

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021



เอกสารไม่ควบคุม

Preventive Maintenance Procedures

FLAME SYSTEM section

☒ Section not applicableElectronic components

- ☒ Review and confirm instrument configuration data in SVD
- ☒ Confirm power supply voltages using the **SVD Power Supply diagnostic**.
- ☒ For Dual Beam instruments - Confirm RBC frequency using the **SVD RBC frequency diagnostic**.

Mechanical components

- ☒ Check the burner adjuster controls for complete and free movement. If the burner adjuster needs lubrication, use Molykote 321 or mineral-based molybdenum disulphide grease.
- ☒ Run SVD tests to exercise all motor drives over the full range of their travel:
 - ☒ Monochromator drive
 - ☒ Slit drive
 - ☒ Lamp selector
 - ☐ ABA *N/A*

Optics components

- ☒ Check that external optical surfaces are clean – Clean or replace as required.
- ☒ Use SVD and perform **Mono Wavelength Correction**.
- ☒ Use SVD and perform **Slit Calibration**.
- ☒ Use SVD and perform **Grating Squareness Diagnostic**.
- ☒ Use SVD and perform **Zero Order Offset/Mono Correction**.
- ☒ Use SVD and perform **Wavelength Repeatability**.
- ☒ Physically inspect selected HC lamps (customer to supply per their choice) and measure the % Gain for each lamp. Advise customer if lamps are showing emission degradation due to age.
- ☒ Check that the signal energy of the D2 and HC lamps track properly. Advise customer if their D2 lamp is showing emission degradation due to age.

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021



เอกสารไม่ควบคุม

Sample Introduction and Atomization

- ☒ Inspect the burner interlock plate to ensure that the interlock pin is secure and correct for the burner type.
- ☒ Clean the burner slot with a clean white card.
- ☒ Check the uniformity of the slot width.
- ☒ Clean the burner if required.
- ☒ Change the burner o-ring.
- ☒ Clean the nebulizer, spray chamber and liquid trap.
- ☒ Change all o-rings and seals in the nebulizer, nebulizer block and spray chamber.
- ☒ Check that the pressure relief bung releases readily.
- ☒ Change o-rings on the fuel and oxidant delivery barbs
- ☒ Leave the liquid trap EMPTY and verify the flame will not ignite in this state.
- ☒ Refill liquid trap and check that overflow drains freely into the drain/waste tube.
- ☒ Check the drain/waste tube for good drainage. It should not have tight bends, kinks or loops and the lower end must be above the liquid level in the waste vessel
- ☒ Check and clean the igniter electrode

Gas handling components and safety interlocks

- ☒ Pressure test for leaks
- ☒ Leak test gasbox internal components and connections
- ☒ Check safety interlock status and operation using the **SVD interlock monitoring diagnostic**.

Analytical performance for Flame systems

- ☒ Ignite a flame.
- ☒ Check that you can adjust the nebulizer uptake rate from 4 to 6.5 mL per minute.
- ☒ Optimize the instrument ready to perform Cu sensitivity test.
- ☒ Create a manual method to perform a Basic Cu ABS test - "Final Performance Testing"
- ☒ Run a PM completed sensitivity test for a 5 ppm copper sample and record the results in the AA PM Performance test results and measurements table.

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021



เอกสารไม่ควบคุม

FURNACE SYSTEM section

☒ Section not applicableElectronic components

- ☐ Review and confirm instrument configuration data in SVD
- ☐ Confirm power supply voltages using the **SVD Power Supply diagnostic**.

Mechanical components

- ☐ Run SVD tests to exercise all motor drives over the full range of their travel:
 - ☐ Monochromator drive
 - ☐ Slit drive
 - ☐ Lamp selector

Optics components


- ☐ Check that external optical surfaces are clean – Clean or replace as required.
- ☐ Use SVD and perform **Mono Wavelength Correction**.
- ☐ Use SVD and perform **Slit Calibration**.
- ☐ Use SVD and perform **Grating Squareness Diagnostic**.
- ☐ Use SVD and perform **Zero Order Offset/Mono Correction**.
- ☐ Use SVD and perform **Wavelength Repeatability**.
- ☐ Physically inspect selected HC lamps (customer to supply per their choice) and measure the % Gain for each lamp. Advise customer if lamps are showing emission degradation due to age.

Gas handling, water system and workhead component checks

- ☐ Inspect the GTA workhead gas hoses and connections for leaks.
- ☐ Pressure test for gas leaks
- ☐ If the cooler system is accessible (stand-alone) check for correct operation and coolant/water level – this includes any temperature and pressure settings plus filter cleaning (air flow and water).
- ☐ Inspect the GTA workhead water hoses and connections for leaks.
- ☐ Check all graphite components and replace if necessary.

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021

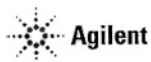

เอกสารไม่ควบคุม

- ☐ Tube
- ☐ Electrodes
- ☐ Shroud

- ☐ Check and clean the end windows on the workhead.
- ☐ Check safety interlock operation.

Analytical performance for Furnace systems

- ☐ Optimize the instrument ready to perform Cu sensitivity test.
- ☐ Run the sensitivity test for a 25 ppb copper sample and record the results in the results table.

PSD autosampler accessory for Furnace systems


- ☒ Section NOT Applicable
- ☐ Check condition of the PSD capillary – replace if necessary
- ☐ Check condition and operation of PSD syringe – ensure it does not have air locks and bubbles.
- ☐ Change PSD rinse bottle o-ring.
- ☐ Check and clean the rinse vessel.
- ☐ Check the drain tube for good drainage. It should not have tight bends, kinks or loops and the lower end must be above the liquid level in the waste vessel.
- ☐ Ensure that the waste vessel is suitable for use with the furnace system.

Sample introduction pump system (SIPS) accessory

- ☒ Section NOT Applicable
- ☐ Re-torque screws securing the hubs, presser arms and pump rotors.
- ☐ Adjust each roller so that it rotates freely.
- ☐ Wipe clean the pump rotor rollers and pump bands with a dry clean cloth.
- ☐ Ensure that the presser arms and the surfaces near the pump are free from dirt and spills.
- ☐ Remove the pump module rear cover and check for the incursion of liquids and any signs of corrosion.
- ☐ Re-torque the nuts that fasten the motor mounting plates to the chassis.
- ☐ Check clips securing the diluents holder and replace if necessary.
- ☐ Disconnect, clean T-piece, and reassemble the tubing using the following steps.

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021

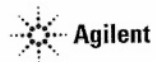

เอกสารไม่ควบคุม

- ☐ Remove the T-piece by disconnecting the pump tubes, the pump bands and all other tubing.
- ☐ Place the T-piece in an ultrasonic bath containing strong detergent 1-5% Decon 30 or similar, for approximately 5-10 minutes.
- ☐ Wash the T-piece under a tap with a strong flow of water.
- ☐ Rinse with distilled water through all of the inlets in the reverse direction to normal sample flow.
- ☐ Reassemble.

Sample preparation system (SPS 4) accessory☒ Section NOT Applicable

The Agilent SPS 4 autosampler is designed to need minimal maintenance.

The following maintenance requirements are suggested to maintain the performance of the autosampler.

- ☐ Cleaning the spill tray, rack location mat, end frames and chassis accessories with a damp soft cloth and diluted mild detergent.
- ☐ Cleaning the autosampler cover panels with domestic window cleaner.
- ☐ Checking the X- axis and Z- axis drive belts for cracks, splits, damaged teeth, excessive fraying, color changes or degradation from fumes..
- ☐ Check the X- axis, Theta- axis and Z- axis FFC cables for cracks, incorrect positioning, damaged edge or damaged connectors.

NOTE: The autosampler requires no extra lubrication throughout its lifetime.


For further details refer to the SPS 4 service manual G8410-90050.

Sample preparation system (SPS 3) accessory☒ Section NOT Applicable

- ☐ Check the x-axis and z-axis timing belts – Replace if there is any cracks, splits or color deterioration and belt tension.
- ☐ Check belt tensions - adjust if required
- ☐ Check the lubrication pad for single x-axis shaft. If pad is dry or customer has observed any vibration or erratic movements of the x-axis carriage, add 1 mL of Dow Corning 200 @ Fluid, 200 CS into the well.
- ☐ Check the auto-sampler ability to find tube positions - Calibrate if required.
- ☐ Clean the exterior surfaces of the accessory with soft lint free cloth. This cloth can be dampened with warm water or a mild detergent. Do not use organic solvents or abrasive cleaning agents.

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021

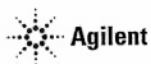

เอกสารไม่ควบคุม

Vapor generation accessory VGA (hydride generator)☒ Section NOT Applicable

- ☐ Inspect VGA gas supply hose.
- ☐ Inspect/replace VGA pump tubing.
- ☐ Check low gas pressure interlock setting – adjust if required.
- ☐ Check precision orifice gas flow setting – adjust if required.
- ☐ Check gas regulator pressure to 46 psi (325 kPa) – adjust if required.
- ☐ Clean the exterior surfaces of the accessory with soft lint free cloth. This cloth can be dampened with warm water or a mild detergent. Do not use organic solvents or abrasive cleaning agents.

UltrAA lamp accessory (external)☒ Section NOT Applicable

- ☐ Check the condition of the power cable.
- ☐ Clean the exterior surfaces of the accessory with soft lint free cloth. This cloth can be dampened with warm water or a mild detergent. Do not use organic solvents or abrasive cleaning agents.

Restore System

- ☒ If you have altered the customer's instrumentation during the course of PM, restore to the original status to allow the customer to conduct their normal activities (e.g., reload the customer's method.)

Guidance

If the PM service is performed prior to a qualification service, then use the qualification procedure as a guide for final instrument set up and checkout.

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021



เอกสารไม่ควบคุม

Signature Page**Service Review**

- ☒ Attach available reports/printouts of all tests to this documentation.
- ☒ Record the Preventive Maintenance service activity in the customer's records/logbook.
- ☒ Update/reset instrument maintenance counters as appropriate.
- ☒ Affix the PM sticker to the system or instrument logbook based on the customer's request.
- ☒ Complete the Service Engineer Comments section if there are additional comments.
- ☒ Review this service, parts replaced, and test results obtained with the customer.
- ☒ If the instrument firmware was updated, record the details of the change in the Service Engineer's Comments box or if necessary, in the customer's IQ records.

Test Results

Test Description		
Flame optics PMT Gain test		
For copper at 324.8 nm, 4 mA, 0.5 nm slit width	< 55 %	44 %
Flame performance test with 5 ppm copper sample		
Air /acetylene, mixing paddle removed	Abs value > 0.5	0.7401 Abs
Air /acetylene, mixing paddle installed, 10 replicates	%RSD < 1.0	0.5 % RSD
Deuterium furnace optics PMT Gain test		
For copper at 324.8 nm, 4 mA, 0.5 nm slit width	< 55 %	N/A
Deuterium furnace performance test with 25 ppb copper sample (324.8 nm)		
Precision %RSD	≤ 4.0%	N/A
Abs value	≥ 0.15	N/A
Zeeman furnace analytical performance: 25 ppb copper sample (327.4 nm)		
Precision %RSD	≤ 4.0%	N/A
Abs value	≥ 0.10	N/A
MSRP	≥ 70 %	N/A

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021



เอกสารไม่ควบคุม

AA consumable and parts list table

Part Description		Product/Model # where used	PM supplied or Consumable	Instrument-Type
Test Solution – Cu 5ppm solution	6610030100	50 55 140 240 280	PM supplied	Common
Test Solution - Blank solution	5190-7001	50 55 140 240 280	PM supplied	Common
Copper, 1000 ug/ml, 100ml	5190-8279	50 55 140 240 280	*	Common
Kit, Mk 7 O-rings, aqueous, complete set	9910093400	50 55 140 240 280	PM supplied	Flame
Organic Kit	9910093500	50 55 140 240 280	PM supplied	Flame
Wire Nebulizer Cleaning	9910024700	50 55 140 240 280	consumable	Flame
Tubing-Capillary Std Nebs	9910024800	50 55 140 240 280	consumable	Flame
Capillary Tube Hivac Neb (3) (organics only)	9910044000	50 55 140 240 280	consumable	Flame
Glass impact beads (5/pk)	9910025700	50 55 140 240 280	consumable	Flame
Teflon impact beads (5/pk): (organics only)	9910053300	50 55 140 240 280	consumable	Flame
Burner cleaning strip (100/pk)	9910053900	50 55 140 240 280	consumable	Flame
Window UV silica – round (right side)	2010082600	50 55 140 240 280	PM supplied	Common
Window UV silica – rectangular (left side)	2010082500	50 55 140 240 280	PM supplied	Common
Pad adhesive window (round)	4910012700	50 55 140 240 280	PM supplied	Common
Pad adhesive window (rectangular)	4910012800	50 55 140 240 280	PM supplied	Common
Electrode kit (1 pr) (D2)	6310003400	GTA120	PM supplied	Furnace
Shroud (D2)	6310003100	GTA120	PM supplied	Furnace
Zeeman electrode kit (1 pr)	6310003500	GTA120	PM supplied	Furnace
Zeeman shroud	6310003600	GTA120	PM supplied	Furnace
O-ring PSD rinse bottle	6910025900	PSD120	PM supplied	Furnace

* For engineers who only service AA instruments 5190-8279 can be used as a cheaper alternative for 6610030100.

Items classified as PM supplied in the above table are included in the standard PM

Those classified as consumable should be provided by the customer or charged to the customer if supplied by the Agilent service engineer.

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021



เอกสารไม่ควบคุม

Service Engineer Comments (optional)

If there are any specific points you wish to note as part of performing the installation or other items of interest for the customer, please write in this box.

Service Completion

Service request number 6006371115 Date service completed 24 January 2024
Agilent signature Worawit T. Customer signature Janida
Total number of pages in this document 13

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021



เอกสารไม่ควบคุม

Power Supply:

Averaging Period: 30.0
Datapoint Count: 20

	Lower Limit (V)	Actual (V)	Upper Limit (V)	Result:
12.00 V Rail	10.80	12.19	13.20	Passed
-12.00 V Rail	-13.20	-11.90	-10.80	Passed
5.00 V Rail	4.50	5.05	5.50	Passed
310.00 V Rail	279.00	320.00	341.00	Passed

เอกสารไม่ควบคุม

SVD Results Report



Report ID: Diagnostic Start Time: 1/24/2024 9:41:24 AM Diagnostic End Time: 1/24/2024 10:10:55 AM

Customer:

Service Engineer: Worawit T.

Address:

Contact Details:

Instrument Configuration

Configuration:

Serial Number: MY13160001 Turret Type: Automatic
Instrument Model: Varian AA140/240/280 Number Of Lamps: 4
Flame Instrument: True Mono Type: Automatic
Furnace Instrument: True Gasbox Type: Y Gas Box
Zeeman Present: False Auto Burner Adjuster: False
Internal Zeeman: False Mains Frequency: 50
Internal UltraAA: False Firmware Version: 2.11
Optics Type: Double Beam Photomultiplier Type: Normal(900nm)
D2 BG Correction Fitted: True PWB Version: 45
Boot Block Version: 1.09

EEPROM Data:

Instrument Run Hours: 62609.832 D2 Run Hours: 49136.000
Zero Wavelength Offset: 30.148 D2 Serial Number: not set !
Mono Correction: 0.765 D2 Install Date: 1/1/1970
Flame Hours: 29802.416 D2 Original Intensity: 1.000
D2 Last Intensity: 475.000

Frequency:

Averaging Period: 30.0
Datapoint Count: 20

Upper Limit: 51.00 Highest Measured Frequency: 50.00
Average Frequency: 50.00
Lower Limit: 49.00 Lowest Measured Frequency: 50.00

Result: **Passed**

Report Generated At: 1/24/2024 10:11:18 AM

1

SVD Results Report



เอกสารไม่ควบคุม

Optics

Beam Balance:

Lamp Type: Copper
Lamp Socket Used: 3

Peak Selected: 324.80
Lamp Alignment: **Performed**

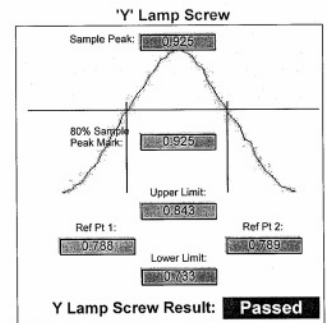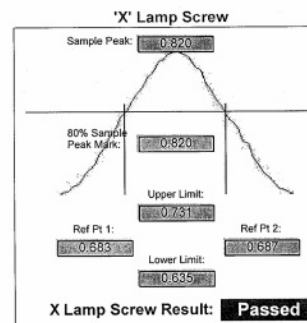

Grating Squareness:

Lamp Element(s): Copper
Lamp Turret Position: 3
Lamp Current(mA): 4.00
Slit Width(nm): 0.5
1st Order Wavelength(nm): 324.80
Lamp Alignment: **Performed**

	Lower Limit (nm)	Actual (nm)	Upper Limit (nm)	Result:
Zero Order	-0.10	0.00	0.10	Passed
First Order	324.45	324.75	325.15	Passed
Second Order	649.23	649.52	649.97	Passed

เอกสารไม่ควบคุม

เอกสารไม่ควบคุม

Wavelength Repeatability:

Lamp Used: Copper

Lamp Current(mA): 4

Peak Used(nm): 324.750

Slit Width(nm): 0.2

Connected to Socket: 3

Slit Height: Normal

Lamp Alignment: Performed

324.888

Upper Limit(nm)

(Approach from Zero Order)

(Approach from end)

Sample 1: 324.828

Sample 2: 324.828

Sample 3: 324.828

Sample 4: 324.823

Sample 5: 324.823

Sample 6: 324.823

Sample 7: 324.823

Sample 8: 324.823

Sample 9: 324.823

Sample 10: 324.823

Mean: 324.825

Standard Deviation: 0.002

Result: Passed

Auto Lamp Recognition:

Lamp 1: Uncoded Lamp/Not Connected

Lamp 5: Not Supported

Lamp 2: 87 - Silver/Cadmium/Lead/Zinc(UltrAA) (Ag/C

Lamp 6: I/Not Supported

Lamp 3: 14 - Copper (Cu)

Lamp 7: Not Supported

Lamp 4: Uncoded Lamp/Not Connected

Lamp 8: Not Supported

Result: Passed

GTA Temperature Monitoring:

Not Performed

Notes:

PM 24 Jan 2024

Signatures:

David

24/1/24

Date

Worawit T.

24/1/24

Date

Mechanical

Wavelength Drive: Passed

Slit Drive: Passed

Turret Drive: Passed

Auto Burner Adjuster Drive: Untested

Miscellaneous

Signal Processing Linearity:

Calculate Mode: New Calc Mode

	Lower Limit	Actual	Upper Limit	Result:
S0	114	261	297	Passed
S1	156	165	191	Passed
S2	271	296	332	Passed
S3	474	507	579	Passed
S4	825	918	1008	Passed
S5	1435	1528	1754	Passed
S6	2498	2769	3053	Passed
S7	4347	4752	5313	Passed

Interlocks:

Burner Fitted: Working

Flame Detect: Working

N2O Burner Fitted: Untested

GCU Active: Working

Flame Shield Closed: Working

Oxidant Pressure: Working

Gas Control Fitted: Untested

Oxidant Changeover: Untested

Pressure Release Bung Fitted: Working

Ignition: Working

Liquid Trap Fitted: Working

Sequential by time report

1/24/2024 11:46 AM

SpectrAA

Page 1 of 1

Analyst

Date Started

Worksheet

Comment

Methods

Computer name

Serial Number:

1/24/2024 11:39 AM

GMT: 1/24/2024 4:39 AM

Cu 5 PPM Sense check

Cu

DESKTOP-R9UFRS

MY13160001

Method: Cu (Flame)

Sample ID	Conc. mg/L	%RSD	Mean Abs
CAL ZERO	0.000	55.0	0.0003
	Readings		
	0.0002	0.0002	0.0004
			1/24/2024
STANDARD 1	5.000	1.7	0.7419
	Readings		
	0.7274	0.7515	0.7468
			1/24/2024

Abs

Linear - Cal. Set 1

0.74

0.60

0.40

0.20

0.00

0.000

2.000

4.000

5.500

Cu mg/L

Curve Fit

Characteristic Conc

r

Calculated Conc

Residuals

= Linear

= 0.028 mg/L

= 1.0000

= 0.000 5.000

= 0.000 0.000


Abs = 0.14833 x C + 0.00026

Sample 001	4.988	0.7	0.7401
	Readings		
	0.7454	0.7399	0.7349
			1/24/2024

Analyst
Date Started 1/24/2024 11:47 AM GMT: 1/24/2024 4:47 AM
Worksheet Cu 5 PPM Precision
Comment
Methods Cu
Computer name DESKTOP-RSUFFS
Serial Number: MY13160001

Method: Cu (Flame)

Sample ID	Exp Abs	%RSD	Mean Abs
Cu Precision	0.723	0.5	0.7232
Readings			
0.7221	0.7195	0.7226	0.7283
0.7201	0.7213	0.7266	0.7174
			1/24/2024

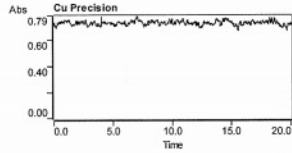

เอกสารไม่ควบคุม

UNITED ANALYST AND ENGINEERING CONSULTANT COMPANY Ltd.

Automatic Mercury Analyzer

Model RA-4500

Preventive Maintenance Report

Serial No. : 17780278

Soft version : Ver 2.0.7

ROM version : Ver 2.0.1

Date : 11 July 2023

PM by : Pradit mayong
(Pradit M.)

Approved by : Pathom S.
(Pathom S.)



Coax Group Corporation Ltd.

1131/62,64,325-331 Nakornchaisri road,

Kwang Thanon Nakornchaisri, Dusit, Bangkok 10300 Thailand

Tel. 02-2435263, 02-6682436 Fax. 02-2437386

เอกสารไม่ควบคุม

Inspection result

ITEM	STANDARD	RESULT	JUDGE
1. Self Check	1.1 Heating	PASS	OK
	1.2 Cooling	PASS	OK
	1.3 Leak	PASS	OK
	1.4 Optical system	PASS	OK
	1.5 Drift	PASS	OK

2. Analytical curve inspection (AREA)			
2.1 No Pretreatment (Low Conc.)	Correlation coefficient	1.0000	OK
	(r) ≥ 0.9990		


3. Repeatability (AREA)			
3.1 No Pretreatment 100ppb, n=5			
	1. 99.12	ppb	
	2. 101.48	ppb	
	3. 101.24	ppb	
	4. 102.34	ppb	
	5. 101.92	ppb	
	C.V. ≤ 5%	1.23%	OK

4. Blank	Below 1.0 (AREA)	0.2062	OK
----------	------------------	--------	----

เอกสารไม่ควบคุม


Title : Preventive Maintenance RA-4500 sn:17780278
Date : 7/11/2023
Name : Coax Group
Memo : Calibration Curve 0-10ng

Calib

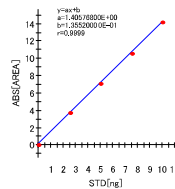

STD

No.	STD [ppb]	SVOL [mL]	CVOL [mL]	DVOL [mL]	STD [ng]	AREA [ON]	MEAS [ng]	Dev [%]	Note
1	100.000	0.000	5.000	5.000	0.000	0.0859	-0.0353	-	
2	100.000	0.025	5.000	5.000	2.500	3.7687	2.5845	3.4	
3	100.000	0.050	5.000	5.000	5.000	7.1028	4.9562	0.9	
4	100.000	0.075	5.000	5.000	7.500	10.6441	7.4753	0.3	
5	100.000	0.100	5.000	5.000	10.000	14.2203	10.0193	0.2	

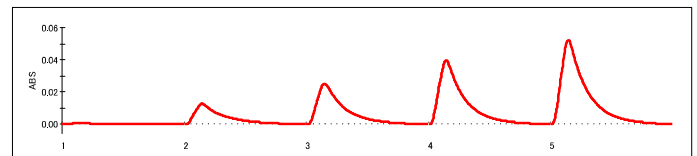

SMP


No.	NAME	SVOL [mL]	CVOL [mL]	DVOL [mL]	AREA [ON]	MEAS [ng]	CONC [ug/L]	Note
1	hg 100 ppb	0.050	5.000	5.000	7.1027	4.9561	99.122	
2	hg 100 ppb	0.050	5.000	5.000	7.2687	5.0742	101.484	
3	hg 100 ppb	0.050	5.000	5.000	7.2514	5.0619	101.238	
4	hg 100 ppb	0.050	5.000	5.000	7.3285	5.1168	102.336	
5	hg 100 ppb	0.050	5.000	5.000	7.2996	5.0962	101.924	

Statistics

No.	NAME	TRY	AV [ug/L]	SD [ug/L]	Cv [%]
1	hg 100 ppb	5	101.2208	1.246264	1.23

เอกสารไม่ควบคุม

NIC NIPPON INSTRUMENTS CORPORATION

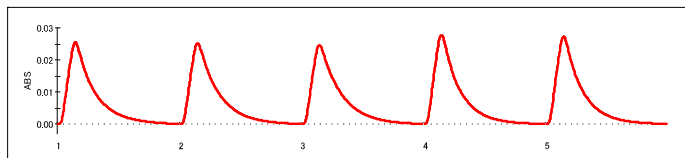


Self Check

Heat check: PASS!! (24.9degC[05:00] -> 28.9degC[02:31])
 Sensor check: PASS!! (78- 18= 60)
 Leak check: PASS!! (0.17L/min)
 Sig/Ref check: PASS!! (Sig: 3.73V, Ref: 3.94V)
 Drift check: PASS!! (-0.0027882 - -0.0032876 = 0.0004993)

Title : Preventive Maintenance RA-4500 sn:17780278
 Date : 7/11/2023
 Name : Coax Group
 Memo : Blank

SMP

No.	NAME	SVOL [mL]	CVOL [mL]	DVOL [mL]	AREA [ON]	MEAS [ng]	CONC [ug/L]	Note
1	Blank				0.2062	0.0503		

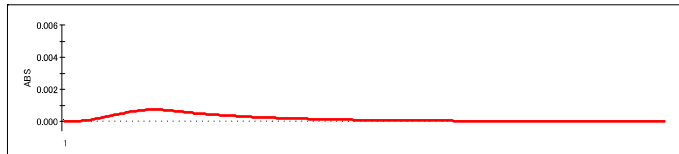

เอกสารไม่ควบคุม

NIC NIPPON INSTRUMENTS CORPORATION

เอกสารไม่ควบคุม

NIC NIPPON INSTRUMENTS CORPORATION

DQE Services Co.,Ltd.
 32 Soi Ladprao-Wanghin 55, Ladprao-Wanghin Rd., Ladprao, Bangkok 10230
 Phone : +66 (0)2 538 2054, Email : dqeservicesinfo@gmail.com

CERTIFICATE OF CALIBRATION

Certificate No. : SP23-021 Page 1 of 5

Customer : United Analyst and Engineering Consultant Co.,Ltd. (Head Office)

Address : 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260

Location of calibration : Laboratory 315

Equipment : UV-Vis Spectrophotometer

Manufacturer : Agilent Technologies

Model : Cary 60

Serial No. : MY15410009

ID No. : N/A

Received Date : 20 May 2023

Calibration Date : 20 May 2023

Issue Date : 23 May 2023

Condition Instrument : Good

Calibrated by : (Mr. Tanasut Rittidach) Technical Manager

Approved by : (Ms. Chonthicha Sangngern) Quality Manager

The calibration result is applied only to the above calibrated item and was found accurate as shown on date and place of calibration only.

The measurement capability of the laboratory and its traceability to recognized national standards and to the unit of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the DQE Services Co., Ltd.

DQE Services Co.,Ltd.
 32 Soi Ladprao-Wanghin 55, Ladprao-Wanghin Rd., Ladprao, Bangkok 10230
 Phone : +66 (0)2 538 2054, Email : dqeservicesinfo@gmail.com

REPORT OF CALIBRATION

Certificate No. : SP23-021 Page 2 of 5

Environment Condition : Ambient Temperature 25 ± 5 °C
 Relative humidity 55 ± 20 %RH

Calibration method : In-house method CP-01 Based on ASTM E275-08

Certified Reference Materials :

Material	Serial No.	Certificate No.	Due date
Absorbance Standard set	25760	95935	22 October 2023
Absorbance Standard set	25757	95929	22 October 2023
Wavelength Standard set	25806	95916	22 October 2023
Wavelength Standard set	25758	95915	22 October 2023

Traceability This certification is traceable to the International System of Unit maintained at National - Institute of Standards and Technology (NIST) through Starna Scientific Limited

Spectral Band Width of UUC : 1.5 nm.

Scan Speed of UUC : 60 nm/min

Scan Interval of UUC : 0.15 nm.

Resolution of UUC : Photometric 0.0001 Abs.

Wavelength 0.1 nm.

เอกสารไม่ควบคุม

เอกสารไม่ควบคุม



REPORT OF CALIBRATION

Certificate No. : SP23-021

Page 3 of 5

Calibration Results : Without adjustment

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor k
420	0.0000	0.0000	0.0000	0.0028	2.00
	0.5787	0.5742	0.0045	0.0031	2.00
	1.0490	1.0423	0.0067	0.0029	2.00
	2.1900	2.1847	0.0053	0.0075	2.00
440	0.0000	0.0000	0.0000	0.0028	2.00
	0.5607	0.5577	0.0030	0.0034	2.00
	1.0247	1.0234	0.0013	0.0035	2.00
	2.1229	2.1171	0.0058	0.0088	2.00
465	0.0000	0.0000	0.0000	0.0028	2.00
	0.5236	0.5184	0.0052	0.0029	2.00
	0.9634	0.9607	0.0027	0.0029	2.00
	1.9763	1.9715	0.0048	0.0081	2.00
546.1	0.0000	-0.0001	0.0001	0.0028	2.00
	0.5191	0.5159	0.0032	0.0031	2.00
	1.0003	0.9980	0.0023	0.0033	2.00
	1.9987	1.9917	0.0070	0.0087	2.00
590	0.0000	0.0000	0.0000	0.0028	2.00
	0.5523	0.5501	0.0022	0.0030	2.00
	1.0809	1.0808	0.0001	0.0030	2.00
	2.0391	2.0336	0.0055	0.0081	2.00
635	0.0000	0.0000	0.0000	0.0028	2.00
	0.5601	0.5585	0.0016	0.0031	2.00
	1.0512	1.0485	0.0027	0.0030	2.00
	1.9294	1.9317	-0.0023	0.0083	2.00

FM-708-02 R01 1/11/2021

เอกสารไม่ควบคุม



REPORT OF CALIBRATION

Certificate No. : SP23-021

Page 4 of 5

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor k
235	0.0000	0.0000	0.0000	0.0050	2.00
	0.7478	0.7436	0.0042	0.0058	2.00
257	0.0000	0.0000	0.0000	0.0050	2.00
	0.8686	0.8648	0.0038	0.0064	2.00
313	0.0000	0.0000	0.0000	0.0050	2.00
	0.2912	0.2908	0.0004	0.0052	2.00
350	0.0000	0.0000	0.0000	0.0050	2.00
	0.6448	0.6398	0.0050	0.0058	2.00

FM-708-02 R01 1/11/2021

เอกสารไม่ควบคุม



REPORT OF CALIBRATION

Certificate No. : SP23-021

Page 5 of 5

Wavelength Accuracy :

CRMs Values (nm.)	UUC Reading (nm.)	Correction (nm.)	Uncertainty (nm.)	Coverage factor k
241.72	242.0	-0.28	0.18	2.00
279.45	279.5	-0.05	0.18	2.00
287.81	287.5	0.31	0.18	2.00
334.06	333.5	0.56	0.18	2.00
360.93	360.3	0.63	0.18	2.00
418.59	418.0	0.59	0.18	2.00
445.94	445.3	0.64	0.18	2.00
453.66	453.0	0.66	0.18	2.00
460.02	459.6	0.42	0.18	2.00
536.59	536.4	0.19	0.18	2.00
637.98	638.3	-0.32	0.18	2.00
431.38	431.0	0.38	0.18	2.00
472.50	472.5	0.00	0.18	2.00
513.47	513.5	-0.03	0.18	2.00
528.88	529.0	-0.12	0.18	2.00
573.17	573.0	0.17	0.18	2.00
585.35	585.0	0.35	0.20	2.00
684.40	684.5	-0.10	0.18	2.00
740.72	741.0	-0.28	0.20	2.00
748.55	748.5	0.05	0.18	2.00
807.03	807.0	0.03	0.18	2.00
879.28	879.5	-0.22	0.18	2.00

Remark : - UUC = Unit Under Calibration

- N/A = Not Available

- The result expanded uncertainty of measurement U is stated as the standard uncertainty of measurement multiplied by the coverage factor k, which for a normal distribution corresponds to a coverage probability of approximately 95%

- * Indicates non TISI accredited

- End of Certificate -

FM-708-02 R01 1/11/2021

เอกสารไม่ควบคุม



CERTIFICATE OF CALIBRATION

Certificate No. : SP24-001

Page 1 of 5

Customer : United Analyst and Engineering Consultant Co.,Ltd. (Head Office)

Address : 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260

Location of calibration : Laboratory 213

Equipment : UV-Vis Spectrophotometer

Manufacturer : Hitachi

Model : U-2900

Serial No. : 21E22-009

ID No. : UAE.WAT.051/2564

Received Date : 4 January 2024

Calibration Date : 4 January 2024

Issue Date : 5 January 2024

Condition Instrument : Good

Calibrated by :

Mr. Tanawit Rittidach
Technical Manager

Approved by :

Ms. Chonthicha Sangngern
Quality Manager

The calibration result is applied only to the above calibrated item and was found accurate as shown on date and place of calibration only.

The measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the DQE Services Co., Ltd.

FM-708-02 R01 1/11/2021

เอกสารไม่ควบคุม



REPORT OF CALIBRATION

Certificate No. : SP24-001

Page 2 of 5

Environment Condition : Ambient Temperature 25 ± 5 °CRelative humidity 55 ± 20 %RH

Calibration method : In-house method CP-01 Based on ASTM E275-08

Certified Reference Materials :

Material	Serial No.	Certificate No.	Due date
Absorbance Standard set	25760	115663	25 October 2025
Absorbance Standard set	25757	115638	25 October 2025
Wavelength Standard set	25806	115657	25 October 2025
Wavelength Standard set	25758	115665	25 October 2025

Traceability : This certification is traceable to the International System of Unit maintained at National -

Institute of Standards and Technology (NIST) through Sarna Scientific Limited

Spectral Band Width of UUC : 1.5 nm.

Scan Speed of UUC : 200 nm/min

Scan Interval of UUC : 0.1 nm.

Resolution of UUC : Photometric 0.001 Abs.

Wavelength 0.1 nm.

เอกสารไม่ควบคุม

FM-708-02 R01 1/11/2021



REPORT OF CALIBRATION

Certificate No. : SP24-001

Page 3 of 5

Calibration Results : Without adjustment

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor k
420	0.0000	0.000	0.0000	0.0028	2.00
	0.5780	0.575	0.0030	0.0031	2.00
	1.0484	1.045	0.0034	0.0029	2.00
	2.1876	2.192	-0.0044	0.0080	2.00
440	0.0000	0.000	0.0000	0.0028	2.00
	0.5595	0.558	0.0015	0.0034	2.00
	1.0239	1.023	0.0009	0.0035	2.00
	2.1230	2.125	-0.0020	0.0079	2.00
465	0.0000	0.000	0.0000	0.0028	2.00
	0.5230	0.520	0.0030	0.0030	2.00
	0.9633	0.961	0.0023	0.0029	2.00
	1.9753	1.975	0.0003	0.0070	2.00
546.1	0.0000	0.000	0.0000	0.0028	2.00
	0.5181	0.516	0.0021	0.0031	2.00
	1.0002	0.997	0.0032	0.0033	2.00
	1.9973	1.993	0.0043	0.0084	2.00
590	0.0000	0.000	0.0000	0.0028	2.00
	0.5517	0.550	0.0017	0.0030	2.00
	1.0803	1.079	0.0013	0.0030	2.00
	2.0373	2.032	0.0053	0.0080	2.00
635	0.0000	0.000	0.0000	0.0028	2.00
	0.5591	0.558	0.0011	0.0031	2.00
	1.0518	1.050	0.0018	0.0030	2.00
	1.9274	1.923	0.0044	0.0079	2.00

เอกสารไม่ควบคุม

FM-708-02 R01 1/11/2021



REPORT OF CALIBRATION

Certificate No. : SP24-001

Page 4 of 5

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor k
235	0.0000	0.000	0.0000	0.0050	2.00
	0.7469	0.743	0.0039	0.0057	2.00
257	0.0000	0.000	0.0000	0.0050	2.00
	0.8674	0.862	0.0054	0.0059	2.00
313	0.0000	0.000	0.0000	0.0050	2.00
	0.2919	0.289	0.0029	0.0051	2.00
350	0.0000	0.000	0.0000	0.0050	2.00
	0.6430	0.641	0.0020	0.0055	2.00

เอกสารไม่ควบคุม

FM-708-02 R01 1/11/2021



REPORT OF CALIBRATION

Certificate No. : SP24-001

Page 5 of 5

Wavelength Accuracy :

CRMs Values (nm.)	UUC Reading (nm.)	Correction (nm.)	Uncertainty (nm.)	Coverage factor k
241.72	241.2	0.52	0.18	2.00
279.45	279.0	0.45	0.18	2.00
287.81	287.4	0.41	0.18	2.00
334.06	333.8	0.26	0.18	2.00
360.93	360.6	0.33	0.18	2.00
418.59	418.4	0.19	0.18	2.00
445.94	445.8	0.14	0.18	2.00
453.66	453.4	0.26	0.18	2.00
460.02	459.8	0.22	0.18	2.00
536.59	536.4	0.19	0.18	2.00
637.98	638.0	-0.02	0.18	2.00
431.38	431.2	0.18	0.18	2.00
472.50	472.5	0.00	0.18	2.00
513.47	513.4	0.07	0.18	2.00
528.88	528.9	-0.02	0.18	2.00
573.17	573.4	-0.23	0.18	2.00
585.35	585.2	0.15	0.20	2.00
684.40	684.4	0.00	0.18	2.00
740.72	741.0	-0.28	0.20	2.00
748.55	748.8	-0.25	0.18	2.00
807.03	807.1	-0.07	0.18	2.00
879.28	879.5	-0.22	0.18	2.00

Remark : - UUC = Unit Under Calibration

- N/A = Not Available

- The result expanded uncertainty of measurement U is stated as the standard uncertainty of measurement multiplied by the coverage factor k,

which for a normal distribution corresponds to a coverage probability of approximately 95%

- * Indicates non TISI accredited


- End of Certificate -

เอกสารไม่ควบคุม

FM-708-02 R01 1/11/2021

DQEServices

DQE Services Co.,Ltd.
32 Soi Ladprao-Wanghin 55, Ladprao-Wanghin Rd., Ladprao, Bangkok 10230
Phone : +66 (0)2 538 2054, Email : dqeservicesinfo@gmail.com


MSC 708-T0 17025
CALIBRATION 0494

CERTIFICATE OF CALIBRATION

Certificate No. : SP24-008Page 1 of 5

Customer : United Analyst and Engineering Consultant Co.,Ltd. (Head Office)

Address : 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260

Location of calibration : Laboratory 315

Equipment : UV-Vis Spectrophotometer

Manufacturer : Hitachi

Model : U-1900

Serial No. : 2021-064


ID No. : UAE.WAS.006/2552


Received Date : 16 January 2024

Calibration Date : 16 January 2024

Issue Date : 19 January 2024

Condition Instrument : Good

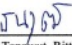Calibrated by :

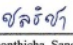(Mr.Tanawut Rittidach)
Technical Manager

Approved by :

(Ms.Chonthicha Sangern)
Quality Manager

The calibration result is applied only to the above calibrated item and was found accurate as shown on date and place of calibration only.


The measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the DQE Services Co., Ltd.

เอกสารไม่ควบคุม

FM-708-02 R01 1/11/2021

DQEServices

DQE Services Co.,Ltd.
32 Soi Ladprao-Wanghin 55, Ladprao-Wanghin Rd., Ladprao, Bangkok 10230
Phone : +66 (0)2 538 2054, Email : dqeservicesinfo@gmail.com


MSC 708-T0 17025
CALIBRATION 0494

REPORT OF CALIBRATION

Certificate No. : SP24-008Page 2 of 5

Environment Condition : Ambient Temperature 25 ± 5 °C

Relative humidity 55 ± 20 %RH

Calibration method : In-house method CP-01 Based on ASTM E275-08

Certified Reference Materials :

Material	Serial No.	Certificate No.	Due date
Absorbance Standard set	25760	115663	25 October 2025
Absorbance Standard set	25757	115638	25 October 2025
Wavelength Standard set	25806	115657	25 October 2025
Wavelength Standard set	25758	115665	25 October 2025

Traceability : This certification is traceable to the International System of Unit maintained at National -
Institute of Standards and Technology (NIST) through Starna Scientific Limited

Spectral Band Width of UUC : 4.0 nm.

Scan Speed of UUC : 200 nm/min

Scan Interval of UUC : 0.1 nm.

Resolution of UUC : Photometric 0.001 Abs.


Wavelength 0.1 nm.

เอกสารไม่ควบคุม

FM-708-02 R01 1/11/2021

DQEServices

DQE Services Co.,Ltd.
32 Soi Ladprao-Wanghin 55, Ladprao-Wanghin Rd., Ladprao, Bangkok 10230
Phone : +66 (0)2 538 2054, Email : dqeservicesinfo@gmail.com


MSC 708-T0 17025
CALIBRATION 0494

REPORT OF CALIBRATION

Certificate No. : SP24-008Page 3 of 5

Calibration Results : Without adjustment

Photometric Accuracy :


Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor k
420	0.0000	0.000	0.0000	0.0028	2.00
	0.5780	0.575	0.0030	0.0031	2.00
	1.0484	1.046	0.0024	0.0029	2.00
	2.1876	2.186	0.0016	0.0080	2.00
440	0.0000	0.000	0.0000	0.0028	2.00
	0.5595	0.558	0.0015	0.0034	2.00
	1.0239	1.024	-0.0001	0.0035	2.00
	2.1230	2.121	0.0020	0.0079	2.00
465	0.0000	0.000	0.0000	0.0028	2.00
	0.5230	0.520	0.0030	0.0030	2.00
	0.9633	0.961	0.0023	0.0029	2.00
	1.9753	1.975	0.0003	0.0070	2.00
546.1	0.0000	0.000	0.0000	0.0028	2.00
	0.5181	0.516	0.0021	0.0031	2.00
	1.0002	0.999	0.0012	0.0033	2.00
	1.9973	1.994	0.0033	0.0084	2.00
590	0.0000	0.000	0.0000	0.0028	2.00
	0.5517	0.550	0.0017	0.0030	2.00
	1.0803	1.080	0.0003	0.0030	2.00
	2.0373	2.032	0.0053	0.0080	2.00
635	0.0000	0.000	0.0000	0.0028	2.00
	0.5591	0.558	0.0011	0.0031	2.00
	1.0518	1.051	0.0008	0.0030	2.00
	1.9274	1.923	0.0044	0.0079	2.00

เอกสารไม่ควบคุม

FM-708-02 R01 1/11/2021

DQEServices

DQE Services Co.,Ltd.
32 Soi Ladprao-Wanghin 55, Ladprao-Wanghin Rd., Ladprao, Bangkok 10230
Phone : +66 (0)2 538 2054, Email : dqeservicesinfo@gmail.com


MSC 708-T0 17025
CALIBRATION 0494

REPORT OF CALIBRATION

Certificate No. : SP24-008Page 4 of 5

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor k
235	0.0000	0.000	0.0000	0.0050	2.00
	0.7469	0.748	-0.0011	0.0057	2.00
257	0.0000	0.000	0.0000	0.0050	2.00
	0.8674	0.865	0.0024	0.0059	2.00
313	0.0000	0.000	0.0000	0.0050	2.00
	0.2919	0.293	-0.0011	0.0051	2.00
350	0.0000	0.000	0.0000	0.0050	2.00
	0.6430	0.641	0.0020	0.0055	2.00

เอกสารไม่ควบคุม

FM-708-02 R01 1/11/2021

DQE Services Co.,Ltd.

32 Soi Ladprao-Wanghin 55, Ladprao-Wanghin Rd., Ladprao, Bangkok 10230

Phone : +66 (0)2 538 2054, Email : dqeservicesinfo@gmail.com

PerkinElmer

For the Better

REPORT OF CALIBRATION

Certificate No. : SP24-008

Page 5 of 5

Wavelength Accuracy :

CRMs Values (nm.)	UUC Reading (nm.)	Correction (nm.)	Uncertainty (nm.)	Coverage factor k
241.54	241.1	0.44	0.18	2.00
279.40	278.9	0.50	0.18	2.00
288.70	288.0	0.70	0.18	2.00
334.22	333.8	0.42	0.18	2.00
361.26	360.8	0.46	0.18	2.00
418.48	418.2	0.28	0.18	2.00
446.70	446.0	0.70	0.18	2.00
453.20	453.1	0.10	0.18	2.00
460.06	459.6	0.46	0.18	2.00
536.90	536.4	0.50	0.18	2.00
637.94	637.6	0.34	0.18	2.00
440.74	440.1	0.64	0.18	2.00
472.22	472.0	0.22	0.18	2.00
513.70	513.5	0.20	0.18	2.00
528.72	528.2	0.52	0.18	2.00
574.60	574.3	0.30	0.18	2.00
585.48	585.0	0.48	0.20	2.00
684.63	684.2	0.43	0.18	2.00
740.27	740.0	0.27	0.20	2.00
748.28	747.8	0.48	0.18	2.00
807.16	806.8	0.36	0.18	2.00
879.70	879.2	0.50	0.18	2.00

Remark : - UUC = Unit Under Calibration

- N/A = Not Available

- The result expanded uncertainty of measurement U is stated as the standard uncertainty of measurement multiplied by the coverage factor k ,

which for a normal distribution corresponds to a coverage probability of approximately 95%

- * Indicates non TISI accredited

- End of Certificate -

เอกสารไม่ควบคุม

PM-708-02 R01 1/11/2021

PinAAcle 900F

Preventive Maintenance Report

Company Name: UNITED ANALYST AND ENGINEERING

Instrument Location: BANGCHAK, PRAKHANONG

BANGKOK, 10260

Instrument Serial No.: PFBS20031902

Date: 26-Jun-2023

เอกสารไม่ควบคุม

PinAAcle 900F Preventive Maintenance (PM)

Company Name:	UNITED ANALYST AND ENGINEERING		
Address (Instrument Location):	BANGCHAK, PRAKHANONG, BANGKOK, 10260		
Serial Number:	PFBS20031902	PM Number:	2/2
Customer Name (if applicable):	K. SATIDA	Telephone Number:	095-558-0049
Customer Support Engineer Name:	K. DUANG	Service Order Number:	WO-02273773
Date PM Performed: (DD-MM-YYYY)	Jun 26, 2023	Next PM Due Date: (DD-MM-YYYY)	Dec 30, 2023
Standard Labor Hours to Complete PM :	5 hours		

Part Number	Release	Publication Date	<div><div>PerkinElmer</div></div>
09370145 Rev.9	A	January 2018	

Component List

Component / Specific Model	Serial #	Configuration Notes

Parts Lists

Parts Included with the PM		
Part Number (if applicable)	Description	Quantity
B0501696	Fan Filters	N/A
N3160156	O-Ring Kits for Sampling Introduction (Stainless Steels Nebulizer)	N/A
N3160157	O-Ring Kits for Sampling Introduction (Plastic Nebulizer)	N/A
N9301714	Replacement Acetylene Filter Cartridge	N/A
TH001022	Replacement Air Filter Cartridge	N/A

Additional Reagents and Standards Required for PM

Part Number (if applicable)	Description	Quality	Batch/Lot #	Expired Date (MM/YY)
N9300183	1000 mg/L Copper Standard	AR	26-87CUV1	30-Jan-2024

Additional Reagents and Standards Required for PM (Customer Support Solution)

Part Number (if applicable)	Description	Quantity	Batch/Lot #	Expiration Date (MM/YY)
N/A	DI Water	250 ml.	AR	AR
N/A	0.5% HNO ₃	250 ml.	AR	AR

เอกสารไม่ควบคุม

เอกสารไม่ควบคุม

Additional Tools Required for PM			
Part Number (if applicable)	Description	Quantity	Serial #
N1013000	0.2A Neutral density filter	1	MGO-252
N1013002	1.0A Neutral density filter	1	MGO-358
03030997	System 2 EDL Driver	1	03030997
N3050605	As System 2 EDL	1	16148
N3050121	Cu Lumina HCL	1	092216-010130
N3050109	Ba Lumina HCL	1	102416-040160
N3050139	K Lumina HCL	1	110716-010060
N3050152	Ni Lumina HCL	1	100516-030190

เอกสารไม่ควบคุม

7. Flame Interlock Check:

Description: Check to ensure that all safety interlocks are closed.

Parameter	Specification	Test Results	Pass/Fail
Flame Sensor	Air/C ₂ H ₂ Flame correctly shuts down	Active	Passed
Drain Sensor	Air/C ₂ H ₂ Flame correctly shuts down	Active	Passed
Nebulizer Sensor	Air/C ₂ H ₂ Flame correctly shuts down	Active	Passed
C ₂ H ₂ Pressure Sensor	Air/C ₂ H ₂ Flame correctly shuts down	Active	Passed
Air Pressure Sensor	Air/C ₂ H ₂ Flame correctly shuts down	Active	Passed
Burner Head Sensor	Choosing Nitrous Oxide as the oxidant should trigger an interlock shuts down	Active	Passed

8. After PM Performance tests:

8.1 Detector Linearity with Barium

Description: Ensures that the detector is linear in the Visible Range.

Parameter	Specification	Certificate Value at 553.6 nm (Abs.)	Test Results	Pass/Fail
1.0 A ND Filter	± 5% from Cert.	0.9798	0.9890	Passed
0.2 A ND Filter	± 5% from Cert.	0.2042	0.1975	Passed

8.2 Baseline Noise at 1.0 Absorbance with Barium

Description: Ensures that a high absorbance will not produce excessive noise.

Parameter	Specification	Results	Pass/Fail
Standard Deviation	≤ 0.010	0.0009	Passed

8.3 AA Baseline Noise with Copper

Description: Check baseline noise.

Parameter	Specification	Results	Pass/Fail
Standard Deviation	≤ 0.001	0.0002	Passed

เอกสารไม่ควบคุม

Procedure Checklist

Use (✓) to check off those steps in the checklist that have been completed.

1. General:

- ✓ Review the instrument performance with the customer and document any recent problems.
- ✓ Inspect the customer log book and make any appropriate PM entries.
- ✓ Perform general inspection of system for cleanliness.

2. PC Instrument Software:

- ✓ Instrument Software user files/databases archived, packed, and/or deleted as needed.

3. Mechanical:

- ✓ Inspect and clean all fans and filters. Replace filters if necessary
- ✓ Inspect all gas lines for leaks and/or wear. Replace if needed.
- ✓ Clean exterior of the instrument.
- ✓ Inspect the burner head, burner chamber, and nebulizer. Clean if needed as stated in the Hardware Guide.
- ✓ Check burner head dimensions with the feeler gauge as stated in the Hardware Guide in the Maintenance chapter section on cleaning the burner head and checking slot width. Replace if out of specification
- ✓ Check the condition of the end cap, burner head, and nebulizer O-rings. Replace if necessary.
- ✓ Check the drain system for signs of wear. Replace worn or damaged parts.
- ✓ Visually check for proper flame conditions when igniting the Air-C₂H₂ and N₂O-C₂H₂ flames (if applicable).

4. Electrical:

- ✓ Inspect PC boards. Clean if necessary.
- ✓ Carefully check all internal and external cable connections.
- ✓ Check instrument firmware revisions upgrade to current levels (if necessary)
- ✓ Run Diagnostics Test within the Advanced function of the Spectrometer page. Check the results in the service log folder in the Spectrometer BM Log Viewer.

5. Optics:

- ✓ Inspect and clean the sample compartment windows, if needed.
- ✓ Inspect optics. Clean or replace if necessary,

6. Gasses:

- ✓ Verify that the Gasses supplied to the instrument are within the pressure and purity specifications found in the PinAAcle 900 Series Pre-installation Checklist SDB.
- ✓ Verify that the acetylene filter and air filter element is dry. Replace if necessary.

เอกสารไม่ควบคุม

8.4 D₂ Background Compensation with Copper

Description: Verifies the instruments ability to compensate for Background absorption.

Parameter	Specification	Results	Pass/Fail
Standard Deviation	≤ 0.010	-0.0062	Passed

8.5 AA-BG Baseline Noise with Copper

Description: Ensures that background correction does not produce excessive noise.

Parameter	Specification	Results	Pass/Fail
Standard Deviation	≤ 0.005	0.0002	Passed

8.6 AA-BG Baseline Noise with Arsenic

Description: Ensures that background correction does not produce excessive noise at a low wavelength.

Parameter	Specification	Results	Pass/Fail
Standard Deviation	≤ 0.005	0.0014	Passed

8.7 Flame Sensitivity

Description: Instrument Sensitivity checked against Copper standard.

Standard Copper Sensitivity	Specification	Results (Abs.)	Pass/Fail
5 mg/L Sensitivity SS Neb (if applicable)	> 0.250 Abs.	NA	Not Applicable
2 mg/L Sensitivity HS Neb (if applicable)	> 0.250 Abs.	0.3467	Passed

10. Review:

- ✓ Review with the customer PM work performed.
- ✓ Review with the customer routine maintenance procedures.
- ✓ Discuss recommended customer supplied materials to have on hand.
- ✓ Attach PM sticker.

เอกสารไม่ควบคุม

Additional Comments

Service Report



Additional Comments Regarding the PM

Review

The preventive maintenance checks and if applicable performance tests for PinAAcle 900F have been completed.

This PinAAcle 900F ☒ Passes ☐ Fails ☐ the preventive maintenance.

Review of Preventive Maintenance:

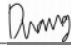Authorized PerkinElmer Representative:		Date: 26-Jun-2023 (DD-MM-YYYY)
Authorized Customer Representative:		Date: 26-Jun-2023 (DD-MM-YYYY)

เอกสารไม่ควบคุม

TO	FOR
Company: United Analyst and Engineering Consultant Co., Ltd. _ Bangkok+HQ Address: 700/2 หมู่ที่ 1 Phrakhanong District, Bangkok, 10260	Work Order Number: WO-00018067 Contact: Karnphong Boonpuang Email: karnphong.b@uaeconsultant.co.th Tel: +66 2763 2828 (7021), +66 8 6347 7390

WORK ORDER INFORMATION			
Top-Level		Order Type	Preventive Maintenance
Installed Product ID	IB-00105024	Billing Type	Paid
Product	SKALAR SAN++ Classic 2SAN59000	PO No.	SSPR2400629
Serial No.	182688	Warranty No.	
		Contract No.	

PRODUCTS SERVICED		
Installed Product Id	Serial Number	Product

PROBLEM DESCRIPTION						
PM 1/1						

Line Number	Engineer	Start Date And Time	End Date And Time	Billable Labor Hour	Billable Travel Hour	Travel KM
WL-00071161	Yongyuth Chanphong	02/20/2024 8:53 AM	02/20/2024 6:07 PM	9.23333		
WL-00092966	Ronnarit Dechnawarat	02/20/2024 8:53 AM	02/20/2024 6:07 PM	9.23333		
Total				18.46666	0	0

! Reach us at DKSH Service-Hotline : +66 2 639 7000
2533 Sukhumvit Road, Bangchak, 10260, Phrakhanong, Bangkok, Thailand
Phone +66 2 639 7000 Fax +66 2 333 1026

Delivering Growth - in Asia and beyond

เอกสารไม่ควบคุม



Line Number	Work Description
WL-00071161	ท่า PM เร็วรอบ
WL-00092966	ท่า PM เร็วรอบ

PARTS CONSUMED		
Part No	Part Description	Quantity

EXPENSES			
Part No	Expense Type	Description	Line Quantity

RECOMMENDED PARTS	
แนะนำอะไหล่ที่ต้องสั่งซื้อเพิ่มเติมดังนี้ , อะไหล่ พารามิเตอร์ Ammonia จำนวน 2 รายการ (9220, 3026) , อะไหล่ พารามิเตอร์ Phenol และ Cyanide จำนวน 6 รายการ (5454, 3028, 3031, 3034, 3036, 3150)	

REMARKS

Travel Time Disclaimer:
Please note that the travel time in this report only includes time taken to reach the installed equipment location. It does not include our engineer's return travel time.

Customer Signature:



Technician: Yongyuth Chanphong
Job Title: Service Manager
Email: yongyuth.yc@dksh.com

Date: 04/10/2024



Job No. WO-00018067

Test Report

Customers	United Analyst and Engineering Consultant Co., Ltd.		
Equipment	Continuous Flow Analyzer	Manufacturer	SKALAR
Controller Mdel	SA5000	Auto Sample Model	SA1052
Controller Serial No.	182688	Auto Sample Serial No.	181729
Date of test	20-Feb-2024	Period	12 Month
Environment temperature	24.7 °C	Humidity	62.2 %RH

Results					
Instrument Checked					
Item	Characteristic	Before	After	Remark	
1	Visual inspect	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail		
2	Power supply (210 - 240 VAC)	<input checked="" type="checkbox"/> 220 VAC	<input checked="" type="checkbox"/> 220 VAC		
3	Computer	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail		
4	Program	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail		
5	Auto sampler	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail		
6	Module holder				
	- Motor pump	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail		
	- Pump tube	<input type="checkbox"/> Pass <input checked="" type="checkbox"/> Fail	<input type="checkbox"/> Pass <input checked="" type="checkbox"/> Fail	*	
	- Air-injection	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail	*	
	- Chemistry manifolds, Switching valve, Coil, Membrane	<input type="checkbox"/> Pass <input checked="" type="checkbox"/> Fail	<input type="checkbox"/> Pass <input checked="" type="checkbox"/> Fail	** , ** *	
7	Detector				
	- Filter	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail		
	- Flow cell	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail		
	- Lamp	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail		
8	Interface	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail		
9	Rinsing valves	<input type="checkbox"/> Pass <input type="checkbox"/> Fail <input type="checkbox"/> N/A	<input type="checkbox"/> Pass <input type="checkbox"/> Fail <input type="checkbox"/> N/A		
10	Temperature / Reactor	<input type="checkbox"/> Pass <input type="checkbox"/> Fail <input type="checkbox"/> N/A	<input type="checkbox"/> Pass <input type="checkbox"/> Fail <input type="checkbox"/> N/A		
11	Flame photometer	<input type="checkbox"/> Pass <input type="checkbox"/> Fail <input type="checkbox"/> N/A	<input type="checkbox"/> Pass <input type="checkbox"/> Fail <input type="checkbox"/> N/A		
12	UPS / Stabilizer	<input type="checkbox"/> Pass <input type="checkbox"/> Fail <input type="checkbox"/> N/A	<input type="checkbox"/> Pass <input type="checkbox"/> Fail <input type="checkbox"/> N/A		

Warning and Error Checked			
Item	Event	Before	After
13	Error list	<input checked="" type="checkbox"/> None <input type="checkbox"/> Appear : _____	<input checked="" type="checkbox"/> None <input type="checkbox"/> Appear : _____

! Reach us at DKSH Service-Hotline : +66 2 639 7000
2533 Sukhumvit Road, Bangchak, 10260, Phrakhanong, Bangkok, Thailand
Phone +66 2 639 7000 Fax +66 2 333 1026

Delivering Growth - in Asia and beyond

เอกสารไม่ควบคุม

DKSH Technology Limited (Head office)
2533 Sukhumvit Road, Bangchak, Phrakhanong, Bangkok, 10260
Phone +66 2 639 7000, Mobile +66 93 813 8681, yongyuth.yc@dksh.com, www.dksh.com

Delivering Growth - in Asia and Beyond



Check with Standard

Item	Characteristic	Before	After	Remark
14	Base Line Test	<input type="checkbox"/> Pass <input type="checkbox"/> Fail <input type="checkbox"/> N/A	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail <input type="checkbox"/> N/A	
15	Detector Signal Test	<input type="checkbox"/> Pass <input type="checkbox"/> Fail <input type="checkbox"/> N/A	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail <input type="checkbox"/> N/A	

Summary of checked

- ☒ The instrument can work normally and efficiently. (เครื่องมือวัดสามารถทำงานได้ปกติและมีประสิทธิภาพ)
- ☐ The instrument can work but it's requiring to maintenance. (เครื่องมือวัดสามารถทำงานได้แต่ต้องบำรุงรักษา)
- ☐ The instrument could not work it's requiring to repair. (เครื่องมือวัดไม่สามารถทำงานได้แต่ต้องการซ่อมบำรุง)

Remark :

* Pump tube และ Air tube เริ่มเสื่อมสภาพ ได้เปลี่ยนอะไหล่ทั้งหมดแล้วตามระยะเวลาใช้งาน
** อะไหล่ 9220 (Manifold T, SiS needle) ขาดไม่สามารถซ่อมได้ และได้เปลี่ยนมาใช้ 5216 ทดแทนชั่วคราว
*** อะไหล่ 5454 (Nipple polyethylene N5) เสีย ได้ทำการเปลี่ยนใหม่แล้ว
หมายเหตุ แนะนำอะไหล่ที่ควรสั่งเพิ่มเติมดังนี้
1. อะไหล่ พารามิเตอร์ Ammonia จำนวน 2 รายการ (9220, 3026)
2. อะไหล่ พารามิเตอร์ Phenol และ Cyanide จำนวน 6 รายการ (5454, 3028, 3031, 3034, 3036, 3150)

Standard Equipment Used

Equipment		Equipment I.D.
Digital multi meter	S/N 57600592	Due date : 8-Jul-2024
Thermo hygrometer	S/N 39520444	Due date : 27-Dec-2024

Test By : 
(Mr. Yongyuth Chanphong)

Approved by : 
(Mr. Eknamong Wankiang)

Position : Supervisor, Technical Service

Position : Manager, Technical Services

DKSH Technology Limited (Head office)
2533 Sukhumvit Road, Bangchak, Phrakhanong, Bangkok, 10260
Phone +66 2 639 7000, Mobile +66 93 813 8681, yongyuth-yc@dksh.com, www.dksh.com

Delivering Growth – in Asia and Beyond

