

ภาคผนวก ข

เอกสารสอบเทียบเครื่องมือวิเคราะห์และตรวจวัดคุณภาพสิ่งแวดล้อม



Certificate No.: T/O 660198

Date of issue : 11-Oct-2023

Equipment Description : Refrigerator
Equipment Model : P1010
Equipment Serial No. : P1010-1020-0005
I.D. No. or Control No. : TNP.LAB.01
Manufacturer : Entech Industrial Solution Co.,Ltd.
Customer Name : TNP ENVIRONMENT CO.,LTD.
Customer Address : 332/173 Moo. 3 Tambon Bang Rak Phatthana, Amphoe Bang Bua Thong,
Nonthaburi 11110
Total pages of certificate : 2 pages
Instrument Receiving Date : 9-Oct-2023
Receiving No. : O-230230
Environmental Conditions : All of the measurement were carried out in the working area
Temperature : (25 ± 15) °C
Humidity : (55 ± 30) % RH
Voltage : (220 ± 22) VAC
Calibration Place : 332/173 Moo. 3 Tambon Bang Rak Phatthana, Amphoe Bang Bua Thong,
Nonthaburi 11110
Calibration Procedure No. : This instrument was calibrated by comparison of reference radiation source standard
according to calibration work instration no WI-CL-18-C

The calibration certificate expended uncertainty of measurement is stated as the standard uncertainty of measurement multiplied by the coverage factor k , which for a normal distribution corresponds to a coverage probability of approximately 95%


*The standard uncertainty of measurement has been determined in accordance with M 3003
The expression uncertainty and confidence in measurement.*

This certificate is applied only to item under test environmental condition.

*This calibration certificate may not be reproduced other than in full except with the permission of the issuing laboratory.
Calibration certificates without signature and seal are not valid.*

This calibration certificate documents are traceability to national standards, which realize the unit of measurement according to the International system of units (SI).

Date of Calibration : 9-Oct-2023



Certificate No. : T/O 660198


The Reference Standard Instrument :-

Instrument	Model	Serial No.	Cert No.	
1) Data logger with RTD Probe	Agilent 34972A	MY41187730	PSL-T 0651-1/66	21-Apr-2024
		MY60008352	PSL-T 0651-3/66	21-Apr-2024

Measured room conditions

Temperature :	Minimum: 30.8 °C	Maximum: 31.9 °C
Humidity :	Minimum: 50.7 %RH	Maximum: 57.2 %RH
Voltage :	Minimum: 219.8 VAC	Maximum: 223.4 VAC
Fresh Air Setting:	off	

Sensor Position :

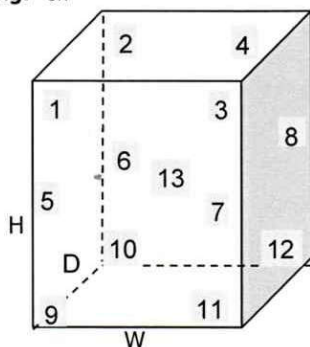

Working Space of chamber :

(Inside Dimensions) W x D x H : 1560 mm x 500 mm x 1380 mm

Sensor Installation Details :

- Sensor Number 1 to 12 installed approximately 50 mm From each wall.
- Sensor Number 13 installed approximately geometric of the chamber.

Results : The measurement results of the calibration were reported in the table below.

(*) Without adjustment

() After adjustment

UUC* Setting	UUC* Reading	Temperature Reading of Standard Sensor								
(°C)	(°C)	Sensor Position								
		1	2	3	4	5	6	7	8	9
		4.02	4.35	4.01	4.20	4.37	4.22	4.17	4.39	4.05
		Sensor Position								
		10	11	12	13					
		4.29	4.30	4.28	4.19					

UUC* Setting	UUC* Reading	Temperature Uniformity	Temperature Stability	Overall Variation	Uncertainty of Measurement	Coverage Factor
(°C)	(°C)	(°C)	(± °C)	(°C)	(± °C)	K
4.0	4.1	1.19	1.08	2.47	1.5	2

UUC* = Unit Under Calibration

Remark :-

- Temperature reading of Standard Sensors shown in the table were taken from the average of Standard reading at each position.
- Temperature Uniformity was calculated from the difference between the maximum and minimum of actual temperature reading from all reference sensors at the same time.
- Temperature Stability was calculated from the maximum stability of nine positions, and formula of Stability is [(Maximum Temperature Value - Minimum Temperature Value) / 2]
- Overall Variation was calculated from the difference between the maximum and minimum measured temperature throughout observation time.

End of Report



CALIBRATION LABORATORY Co., LTD.

2/10-11,14, 55 Soi Prasert Manukit 29 Yaek 4, Prasert Manukit Rd., Ladphrao, Bangkok 10230
Tel. 02-578-0353-4 Fax: 02-578-2672 www.cal-laboratory.com E-mail:sale@cal-laboratory.com



CERTIFICATE OF CALIBRATION

FOR

NOMENCLATURE : pH METER
MANUFACTURER : HORIBA
MODEL / TYPE : LAQUA-PH1100/9615S
SERIAL NO. : B80A0042/9X0B0575
CLID. NO. : 272001452
JOB CONTROL NO. : 230911100397

CUSTOMER : TNP ENVIRONMENT CO., LTD.
332/173 MOO 3 TAMBON BANG RAK PHATTANA,
AMPHOE BANG BUA THONG, NONTHABURI 11110

DATE OF RECEIVED : 11 September 2023

DATE OF ISSUED : 14 September 2023

Report of calibration screening must not be taken in part. Except complete. Without the approval of the Calibration Laboratory Co., Ltd.

REPORT OF CALIBRATION

FOR

NOMENCLATURE : pH METER
MANUFACTURER : HORIBA
MODEL / TYPE : LAQUA-PH1100/9615S
SERIAL NO. : B80A0042/9X0B0575
DATE OF CALIBRATION : 12 September 2023

ENVIRONMENT CONDITIONS :

Temperature : $(25 \pm 2.5) ^\circ\text{C}$

Relative Humidity : $(50 \pm 15) \% \text{ RH}$

PROCEDURE USED :

This instrument was calibrated under procedure No. **CLC-CPCH-01**. The calibration was performed by direct measurement with Certified Reference Material (CRM).

REFERENCE STANDARD USED :

1. pH Standard Solution, NIMT TRM CODE TRM-S-2003, TRM CODE TRM-S-2007.
2. pH Standard Solution, Control Company Catalog Number 06664263,11784256, Lot Number CC752722.

TRACEABILITY :

1. The measurements are traceable to International System of Units (SI) , through National Institute of Metrology (Thailand).
Lot Number. 040822 , 230822. Due Date 26 April 2024.
2. The measurements are traceable to International System of Units (SI) , through Control Company.
Certificate No. 4288-13355261 , Due Date 06 May 2024.

UNCERTAINTY :

The reported expanded uncertainty of measurement is stated as the standard uncertainty of measurement multiplied by the coverage factor complies with the table which for a normal distribution corresponds to a coverage probability of approximately 95 %.

It has been evaluated according to the "Evaluation of the Uncertainty of Measurement in Calibration (EA-4/02 M:2022)"

Certificate No. Q23100397

F3-011-04/01-12

page 2 of 3

CONDITION OF CALIBRATION ITEM : GOOD

MEASUREMENT RESULTS : (X) without adjustment () adjustment

The table in the following gives the calibration results and associated measurement uncertainties of pH meter.

CALIBRATION DATA

pH METER RESULT @ 25 °C

Standard pH Buffer Solution (pH)	pH Meter Reading (pH)	pH Meter Reading (mV)	Correction (pH)	Uncertainty of pH Measurement (\pm pH)	k Factor
4.003	4.01	150.2	-0.007	0.010	2,00
7.000	7.00	-26.1	0.000	0.015	2,06
10.003	10.01	-187.1	-0.007	0.016	2,05

Technical Note. Setting function CAL 3 point (4,7,10).

The Scope of Accredited ANAB Certificate No. ACDM-2814 Version 008 Page 2,3 of 54

This report is valid for the above stated instrument/s only.

End of Certificate

Certificate No. Q23100397

F3-011-04/01-12

page 3 of 3



@clccalibration

CERTIFICATE OF CALIBRATION

FOR

NOMENCLATURE : DIGITAL THERMO-HYGROMETER
MANUFACTURER : EXTECH INSTRUMENTS
MODEL / TYPE : 445814
SERIAL NO. : PONPE5816745
CLID. NO. : 232303263
JOB CONTROL NO. : 230911100396

CUSTOMER : TNP ENVIRONMENT CO., LTD.
332/173 MOO 3 TAMBON BANG RAK PHATTANA,
AMPHOE BANG BUA THONG, NONTABURI 11110

DATE OF RECEIVED : 11 September 2023

DATE OF ISSUED : 15 September 2023

Report of calibration screening must not be taken in part. Except complete. Without the approval of the Calibration Laboratory Co., Ltd.

Certificate No. Q23100396

F3-011-04/01-12

page 1 of 3



@clccalibration

REPORT OF CALIBRATION

FOR

NOMENCLATURE : **DIGITAL THERMO-HYGROMETER**
MANUFACTURER : **EXTECH INSTRUMENTS**
MODEL / TYPE : **445814**
SERIAL NO. : **PONPE5816745**
DATE OF CALIBRATION : **13 September 2023**

ENVIRONMENT CONDITIONS :

Temperature : $(23 \pm 2) ^\circ\text{C}$

Relative Humidity : $(55 \pm 10) \% \text{RH}$

PROCEDURE USED :

This instrument was calibrated under procedure No. **CLC-CPTH-11**. The calibration was performed by using Chilled Mirror Hygrometer and Temperature & Humidity Chamber which maintained by the Calibration Laboratory Co., Ltd.

REFERENCE STANDARD USED :

Chilled Mirror Hygrometer, Edgetech Model Dew Master S/N. 36151.

Temperature & Humidity Chamber, PGC Model 9141-5114 S/N.0802282.

TRACEABILITY :

The measurements are traceable to International System of Units (SI) , through Thunder Scientific Corporation.

Certificate No. 21028, Due Date 09 December 2023.

UNCERTAINTY :

The reported expanded uncertainty of measurement is stated as the standard uncertainty of measurement multiplied by the coverage factor $k = 2,00$ which for a normal distribution corresponds to a coverage probability of approximately 95 %.

It has been evaluated according to the "Evaluation of the Uncertainty of Measurement in Calibration (EA-4/02 M:2022)"

Certificate No. Q23100396

F3-011-04/01-12

page 2

CONDITION OF CALIBRATION ITEM : GOOD

MEASUREMENT RESULTS : (X) without adjustment () adjustment

The table in the following gives the calibration results and associated measurement uncertainties of the measuring digital thermo-hygrometer.

CALIBRATION DATA

1. CORRECTION OF TEMPERATURE

Test point (° C)	Actual Temperature (° C)	DUC Reading (° C)	Correction (° C)	Uncertainty \pm (° C)
20.0	20.01	19.9	+0.11	0.27
25.0	25.01	25.2	-0.19	

2. CORRECTION OF HUMIDITY

STD Temperature (° C)	STD Reading (%RH)	DUC Reading (%RH)	Correction (%RH)	Uncertainty \pm (%RH)
25	50.0	47	+3.0	0.8

Note. The Scope of Accredited ANAB Certificate No. ACDM-2814 Version 008 Page 49 of 54

This report is valid for the above stated instrument/s only.

End of Certificate

Certificate No. Q23100396

F3-011-04/01-12

page 3 of 3



@clccalibration

Certificate of Calibration

Certificate No. : 67-400049-1

Page : 1 of 2

Submitted by : TNP Environment Co., Ltd.

332/173 Moo 3 Bang Rak Phatthana, Bang Bua Thong, Nonthaburi 11110

Equipment : Liquid in Glass Thermometer

Manufacturer : SK

Model : N/A

Range : 0 °C to 200 °C

Resolution : 1 °C

Serial No. : N/A

Immersion : Total

ID No. : TNP.LAB.12

Environment : Ambient Temperature : (23 ± 2) °C

Relative Humidity : (50 ± 15) %

Line Voltage : (220 ± 22) VAC

Date of Received : 26 January 2024

Date of Calibration : 01 February to 02 February 2024

Date of Issue : 02 February 2024

Calibrated by :

Calibration Method : This instrument was calibrated by In-house method comparison technique CAL-M4001 based on ASTM E77-07 by compared with PRT in the liquid bath at the constant controlled temperature.

The temperature scale used was based on ITS-90

Reference Standard Instruments : This certification is traceable to the International System of Units

1. Platinum Resistance Thermometer (PRT)

ID No.	Cert. No.	Due Date	Traceability
400001	TT-0016-22	07 Feb 2024	National Institute of Metrology Thailand (NIMT)

2. Standard Digital Thermometer

ID No.	Cert. No.	Due Date	Traceability
400003	23E1866	01 Jun 2025	National Institute of Metrology Thailand (NIMT)
400004	23E1866	01 Jun 2025	National Institute of Metrology Thailand (NIMT)



Certificate of Calibration

Certificate No. : 67-400049-1

Page : 2 of 2

Result of Calibration : Without Adjustment

UUC Condition As-Received : Good

Function : Temperature measurement

Ice point check : UUC* reading 0 ° C Standard reading 0.8789 ° C

Standard Reading (° C)	UUC Reading (° C)	Correction (° C)	Uncertainty (± ° C)
21.2064	20	1.2	0.31
31.3084	30	1.3	0.31

Remark

UUC : Unit Under Calibration

This result of calibration was found accurate as shown on date and place of calibration only.

This reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%

- oOo -

Certificate no: H/T 670338

Date of issue : 21-Mar-24

Instrument description	:	Thermo-Hygrometer
Instrument model	:	Extech 445815
Instrument serial no.	:	PONPE5899554
ID no. or control no.	:	TNP.LAB.21
Manufacturer	:	Extech Instruments
Probe description	:	-
Probe model	:	-
Probe serial	:	-
Customer name	:	TNP ENVIRONMENT CO.,LTD.
Customer address	:	332/173 Moo 3 Tambon Bang Rak Phatthana, Amphoe Bang Bua Thong, Nonthaburi 11110
 Total pages of certificate	:	 2 Pages
Receiving no.	:	L-241004-1
Receiving date.	:	08-Mar-24
Parameter of calibration	:	Temperature Calibration
Condition of UUC.	:	Used
Ambient condition	:	All of the Measurement were carried out the stabilized laboratory
		Temperature : 23 ± 5 °C
		Humidity : 55 ± 15 %RH
Calibration place	:	17/121 Soi Ngamwongwan 47 Yaek 48, Toongsonghong, Laksi, Bangkok 10210
 Calibration procedure no.	:	 This instrument was calibrated by comparison of indication with the Standard Thermo- hygrometer according to calibration Work Instruction no .WI-CL-11-C

*The calibration certificate expanded uncertainty of measurement is stated as the standard uncertainty of measurent
Multiplied by coverage factor k=2, which for a normal distribution corresponds to a coverage probability of approximately 95%.*

This certificate is applied only to item under test Environmental condition.

*This Calibration Certificate may not be reporduced other than in full except with the permission of the issuing laboratory.
Calibration certificates without signature and seal not valid.*

*This calibration certificate documents are tracebility to national standards, which realize measurement according to the
International System of Units (SI).*

Certificate no: H/T 670338

Standard references

Standard	Reference No.	Vendor	Due Date
ARALAB 300ECP,Fitoclima	S2023070040-001	MIT	07-Jul-24
Thermo HygroPalm HP 23-A	SG-H-00579/66	Success Gateway	16-Aug-24

Measured room conditions

Temperature : 22.1 °C

Humidity : 55.9 %RH

Pressure : 1019.3 mbar

Calibration results (Without Adjustment)

Reference temperature : - °C

Parameter of standard	Standard values	Mean of UUC.	Error	Uncertainty (±)
Temperature (°C)	19.97	20.1	0.13	0.50
Temperature (°C)	25.02	25.2	0.18	0.50
Temperature (°C)	29.99	30.2	0.21	0.50

Remark : -

End of Report

Certificate no: H 670285

Date of issue : 21-Mar-24

Instrument description : Thermo-Hygrometer
Instrument model : Extech 445815
Instrument serial no. : PONPE5899554
ID no. or control no. : TNP.LAB.21
Manufacturer : Extech Instruments
Probe description : -
Probe model : -
Probe serial : -
Customer name : TNP ENVIRONMENT CO.,LTD.
Customer address : 332/173 Moo 3 Tambon Bang Rak Phatthana, Amphoe Bang Bua Thong, Nonthaburi 11110

Total pages of certificate : 2 Pages
Receiving no. : L-241004
Receiving date. : 08-Mar-24
Parameter of calibration : Humidity Calibration
Condition of UUC. : Used
Ambient condition : All of the Measurement were carried out the stabilized laboratory
 Temperature : 23 ± 5 °C
 Humidity : 55 ± 15 %RH
Calibration place : 17/121 Soi Ngamwongwan 47 Yaek 48, Toongsonghong, Laksi, Bangkok 10210
Calibration procedure no. : This instrument was calibrated by comparison of indication with the Standard Thermo- hygrometer according to calibration Work Instruction no .WI-CL-11-C

The calibration certificate expanded uncertainty of measurement is stated as the standard uncertainty of measurent Multiplied by coverage factor $k=2$, which for a normal distribution corresponds to a coverage probability of approximately 95%.

This certificate is applied only to item under test Environmental condition.

This Calibration Certificate may not be reporduced other than in full except with the permission of the issuing laboratory. Calibration certificates without signature and seal not valid.

This calibration certificate documents are tracebility to national standards, which realize measurement according to the International System of Units (SI).

Certificate no: H 670285

Standard refereces

Standard	Reference No.	Vendor	Due Date
Thermo HygroPalm HP 23-A	SG-H-00579/66	Success Gateway	16-Aug-24
Hygrogen 2-XL	Performance	Entech	24-Mar-24

Measured room conditions

Temperature : 22.7 °C **Humidity :** 56.7 %RH **Pressure :** 1013.3 mbar.

Calibration results (Without Adjustment)

Reference temperature : 25.0 °C

Parameter of standard	Standard values	Mean of UUC.	Error	Uncertainty (±)
Humidity (%RH)	35.09	31	-4.09	1.3
Humidity (%RH)	50.03	48	-2.03	1.5
Humidity (%RH)	64.97	68	3.03	1.5

Remark : -

End of Report



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG, BANGKOK 10250
TEL. 0-2717-3000-24 FAX. 0-2719-9484



Certificate of Calibration

Certificate No. : 23M455

Page : 1 of 2

Equipment : Standard Weight Set

Manufacturer: -

Model : Class:F1

Serial No.: 15022021-01

ID No.: TNP.LAB.25

Condition As-Received: Used Item

Received Date: 02 March 2023

Calibration Date: 04 March 2023

Reference: 2303-0104WN

Submitted by: TNP ENVIRONMENT CO.,LTD.

Ambient Temperature: (23 ± 2) °C

Relative Humidity: (50 ± 15) %

Atmospheric Pressure: 1015.25 mbar

This certificate may not be reproduced other than in full,
except with the prior written approval of the head of
Corporate Services 3: Equipment Calibration and Testing Services.

332/173 Moo 3, Bang Rak Phatthana, Bang Bua Thong,
Nonthaburi 11110

Procedure used: Calibration were conducted using in-house calibration procedure CP-M01 according to comparison method against standard weights on the basis of weighings at an average air density of 1.2 kg/m³ and a temperature of 23.4 °C material density of weight is 8000 kg/m³.

Condition of this result of calibration

1.Reference standards instruments :

<u>Instrument</u>	<u>Model</u>	<u>Serial No.</u>	<u>Certificate No.</u>	<u>Due Date</u>
1) Standard Weight Set (E2)	73336	20026	MM-0018-22	28 Feb 2024
2) Standard Weight Set (E2)	73338	20028	MM-0019-22	28 Feb 2024

2.This certificate is not certified for any commercial transaction.

3.The certificate is valid only to the item calibrated on date and place of calibration.

4.This Certification is traceable to the International System of Unit maintained at:-

-National Institute of Metrology Thailand (NIMT)

B 0309848



Cert No.: 23M455

Page: 2 of 2

Result of calibration

Nominal Value	Conventional mass		Uncertainty of Measurement (\pm)	Maximum Permissible error (\pm)
	Before Adjustment	After Adjustment		
200 g	199.99986 g	-	0.30 mg	1.0 mg
100 g	100.00015 g	-	0.16 mg	0.50 mg
50 g	50.00015 g	-	0.10 mg	0.30 mg
20 g	20.000116 g	-	0.080 mg	0.25 mg
10 g	10.000041 g	-	0.060 mg	0.20 mg
5 g	5.000010 g	-	0.050 mg	0.16 mg
2 g	1.999936 g	-	0.040 mg	0.12 mg
1 g	0.999973 g	-	0.030 mg	0.10 mg
200 mg	200.059 mg	200.007 mg	0.020 mg	0.060 mg
100 mg	100.037 mg	99.981 mg	0.016 mg	0.050 mg

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95 %.

-o0o-



a 1151188



THAI CALIBRATION SERVICES CO., LTD.

19/8 Moo 9 Soi Raiking 30 Puttamonthon 5 Rd., Sampran, Nakornpatom 73210

Tel. 0-3439-7682-5 Fax: 0-3439-7687

www.thaical.com E-mail : sale@thaicalibration.com, lab@thaicalibration.com



CALIBRATION CERTIFICATE

Certificate No.S2306518S

page 1 of 2

Customer : TNP ENVIRONMENT CO., LTD.
332/173 Moo 3 Tambon Bang Rak Phatthana,
Amphoe Bang Bua Thong, Nonthaburi 11110

Equipment : Non-automatic weighing instrument (Electronic instrument)

Manufacturer : Shimadzu **Order No. :** 66S2523-1

Model : AP225WD **Ambient temperature :** $(26.9 \pm 5.0) ^\circ\text{C}$

Accuracy class : - **Relative humidity :** $(52.0 \pm 10.0) \%$

Capacity : 10 g / 220 g **Received date :** 21-Jun-2023

Resolution : 0.00001 g / 0.0001 g **Date of calibration :** 21-Jun-2023

Serial No. : D316301848 **Date of issue :** 24-Jun-2023

ID No. : TNP.LAB.30 **Condition of the balance :** Good working conditions

Place of calibration : ห้อง LAB

Calibration method

This instrument was calibrated according to the EURAMET Calibration Guide No. 18.

Condition of reference standard weight

<u>Instrument</u>	<u>Nominal value</u>	<u>Serial No.</u>	<u>Certificate No.</u>	<u>Due-date</u>	<u>Density (kg/m³)</u>
1 Standard weight set	1 mg to 2 kg	15885+15849	M2210001S	8-Oct-2023	7950

Traceability of the reference standard weight

This certificate is traceable to SI unit through Mass Calibration Laboratory Thai Calibration Services Co., Ltd., NSC-ONSC accredited no. Calibration 0189.



THAI CALIBRATION SERVICES CO., LTD.

19/8 Moo 9 Soi Raking 30 Puttamonthon 5 Rd., Sampran, Nakornpatom 73210

Tel. 0-3439-7682-5 Fax: 0-3439-7687

www.thaical.com E-mail : sale@thaicalibration.com, lab@thaicalibration.com



CALIBRATION CERTIFICATE

Certificate No.S2306518S


page 2 of 2

The repeatability of indication

Nominal Value (g)	Standard Deviation of reading (g)	Maximum difference between suscensive reading (g)	n
100	0.000005	0.00001	5
200	0.00005	0.0001	5

The effect of eccentric application of a load on the indication (test load : 100 g)

Position	Balance Reading (g)
Point 1	100.0000
Point 2	100.0002
Point 3	100.0001
Point 4	100.0000
Point 5	100.0001
Eccentric Value	0.0002

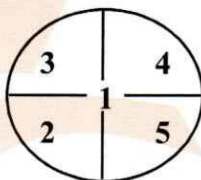

The error of indication

Nominal Value (g)	Value of Reference Standard Weight (g)	Balance Reading (g)	Correction (g)	Uncertainty (±) (g)	k
Unload	0.00000	0.00000	0.00000	0.000016	2.32
0.1	0.10000	0.10003	-0.00003	0.000019	2.10
0.5	0.50000	0.50001	-0.00001	0.000023	2.04
1	1.00001	1.00000	+0.00001	0.000026	2.00
5	5.00000	5.00001	-0.00001	0.000038	2.00
10	9.99999	10.00001	-0.00002	0.000046	2.00
20	20.0000	20.0000	0.0000	0.000085	2.00
50	50.0000	50.0001	-0.0001	0.00011	2.00
100	100.0000	100.0000	0.0000	0.00018	2.00
200	200.0000	200.0004	-0.0004	0.00034	2.00

Remark : Adjustment, External weight nominal value 100 g, Standard weight of Lab

Uncertainty of measurement

The reported expanded uncertainty of measurement is stated as the standard uncertainty of measurement multiplied by the coverage factor (k), which for a normal distribution corresponds to a coverage probability of approximately 95% (confidence level).

This report will certify of the calibrated equipment only.

--End--



THAI CALIBRATION SERVICES CO., LTD.

19/8 Moo 9 Soi Raiking 30 Puttamonton 5 Rd., Sampran, Nakornpatom 73210

Tel. 0-3439-7682-5 Fax: 0-3439-7687

www.thaical.com E-mail : sale@thaicalibration.com, lab@thaicalibration.com



CALIBRATION CERTIFICATE

Certificate No.S2306519S

page 1 of 2

Customer : TNP ENVIRONMENT CO., LTD.
332/173 Moo 3 Tambon Bang Rak Phatthana,
Amphoe Bang Bua Thong, Nonthaburi 11110

Equipment : Non-automatic weighing instrument (Electronic instrument)

Manufacturer : Sartorius **Order No. :** 66S2523-2

Model : SECURA224-1S **Ambient temperature :** $(26.8 \pm 5.0) ^\circ\text{C}$

Accuracy class : - **Relative humidity :** $(52.0 \pm 10.0) \%$

Capacity : 220 g **Received date :** 21-Jun-2023

Resolution : 0.0001 g **Date of calibration :** 21-Jun-2023

Serial No. : 0041305301 **Date of issue :** 24-Jun-2023

ID No. : TNP.LAB.31 **Condition of the balance :** Good working conditions

Place of calibration : ห้อง LAB

Calibration method

This instrument was calibrated according to the EURAMET Calibration Guide No. 18.

Condition of reference standard weight

<u>Instrument</u>	<u>Nominal value</u>	<u>Serial No.</u>	<u>Certificate No.</u>	<u>Due-date</u>	<u>Density (kg/m³)</u>
1 Standard weight set	1 mg to 2 kg	15885+15849	M2210001S	8-Oct-2023	7950

Traceability of the reference standard weight

This certificate is traceable to SI unit through Mass Calibration Laboratory Thai Calibration Services Co., Ltd., NSC-ONSC accredited no. Calibration 0189.

No. 17068



THAI CALIBRATION SERVICES CO., LTD.

19/8 Moo 9 Soi Raiking 30 Puttamonthon 5 Rd., Sampran, Nakornpatom 73210

Tel. 0-3439-7682-5 Fax: 0-3439-7687

www.thaical.com E-mail : sale@thaicalibration.com, lab@thaicalibration.com



CALIBRATION CERTIFICATE

Certificate No.S2306519S


page 2 of 2

The repeatability of indication

Nominal Value (g)	Standard Deviation of reading (g)	Maximum difference between successive reading (g)	n
200	0.00000	0.0000	5

The effect of eccentric application of a load on the indication (test load : 100 g)

Position	Balance Reading (g)
Point 1	100.0000
Point 2	100.0000
Point 3	100.0000
Point 4	99.9998
Point 5	99.9998
Eccentric Value	0.0002

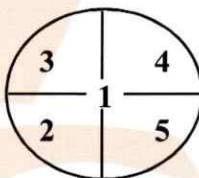

The error of indication

Nominal Value (g)	Value of Reference Standard Weight (g)	Balance Reading (g)	Correction (g)	Uncertainty (±) (g)	k
Unload	0.0000	0.0000	0.0000	0.000082	2.00
0.1	0.1000	0.1000	0.0000	0.000083	2.00
0.5	0.5000	0.5000	0.0000	0.000084	2.00
1	1.0000	0.9999	+0.0001	0.000085	2.00
5	5.0000	5.0000	0.0000	0.000090	2.00
10	10.0000	10.0000	0.0000	0.000094	2.00
20	20.0000	20.0001	-0.0001	0.00011	2.00
50	50.0000	50.0001	-0.0001	0.00013	2.00
100	100.0000	100.0000	0.0000	0.00019	2.00
200	200.0000	199.9998	+0.0002	0.00033	2.00

Remark : Without adjustment

Uncertainty of measurement

The reported expanded uncertainty of measurement is stated as the standard uncertainty of measurement multiplied by the coverage factor (k), which for a normal distribution corresponds to a coverage probability of approximately 95% (confidence level).

This report will certify of the calibrated equipment only.

--End--

CERT.No.: HS-U039F

Harikul Science Co.,Ltd.
 694 Soi Ratchadanivet 24, Pracharatbamphen,
 Samsaennok, Huaikhwang, Bangkok 10310
 Tel: 0-2274-2456 Fax: 0-2274-2443
 Email: info@harikul.com www.harikul.com

Certificate of Calibration

Calibration Date : 20 Jun 23
 Submitted by : TNP ENVIRONMENT COMPANY LIMITED.
 332/173 Moo. 3, Tambon Bang Rak Phatthana,
 Amphoe Bang Bua Thong, Nonthaburi 11110

Model : YSI 4010-2W
 S/N : 22051520
 Probe : YSI 4100
 S/N : 22C102711
 ID NO. : -
 Air Temp ref : S/N. E00522
 Barometric ref : S/N. E00522
 Water Temp ref : S/N. 11431
 Technician : Kittipong M.

Avg Room Temp : 20 °C
 Avg Water Temp : 20 °C
 Air Pressure : 757.00 mmHg
 Salinity : 0 ppt

Calibration Details

Calibration Point	100% air sat. (@20 °C, DO = 9.09 mg/l)	(status)	(status)
Measurement 1 (mg/l)	9.05	(PASS)	-
Measurement 2 (mg/l)	9.05	(PASS)	-
Measurement 3 (mg/l)	9.04	(PASS)	-
Measurement 4 (mg/l)	9.03	(PASS)	-
Measurement 5 (mg/l)	9.04	(PASS)	-
Measurement 6 (mg/l)	9.04	(PASS)	-
Measurement 7 (mg/l)	9.04	(PASS)	-
Measurement 8 (mg/l)	9.03	(PASS)	-
Measurement 9 (mg/l)	9.03	(PASS)	-
Measurement 10 (mg/l)	9.03	(PASS)	-
Mean Measurement	9.04	mg/l	-
Inaccuracy	0.05	mg/l	-
Overall Status	(PASS)		

Manufacturer Specification

Accuracy = +/- 0.2 mg/l

Certificate of Calibration

Certificate No. : 66-400687-1

Page : 1 of 2

Submitted by : TNP Environment Co., Ltd.
332/173 Moo 3 Bang Rak Phatthana, Bang Bua Thong, Nonthaburi 11110

Equipment : Air Chamber (Oven)
Manufacturer : Memmert Model : UF75
Range : N/A °C Resolution : 0.1 °C
Serial No. : B320.0251 ID No. : N/A

Environment : On site calibration was carried out at the Laboratory, TNP Environment Co., Ltd.
Ambient Temperature : (27.0 to 28.0) °C
Relative Humidity : (45 to 50) %
Line Voltage : (228.0 to 230.0) V

Date of Received : 11 December 2023

Date of Calibration : 11 December 2023

Date of Issue : 14 December 2023

Calibrated by : 

Calibration Method : CAL-M4004, TLAS G-20

The temperature scale used was based on ITS-90

Reference Standard Instruments : This certification is traceable to the International System of Units
Standard Digital Thermometer with Thermocouple probe

ID No.	Cert. No.	Due Date	Traceability
400046 & 400023	66-400547-1	04 Apr 2024	National Institute of Metrology Thailand (NIMT)

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full except with the prior written approval of the Calibratech Co.,Ltd.



Certificate of Calibration

Certificate No. : 66-400687-1


Page : 2 of 2

Result of Calibration : Without Adjustment

UUC Condition As-Received : Good

Function : Temperature measurement

This instrument was setting air ventilation at position 0 (close)

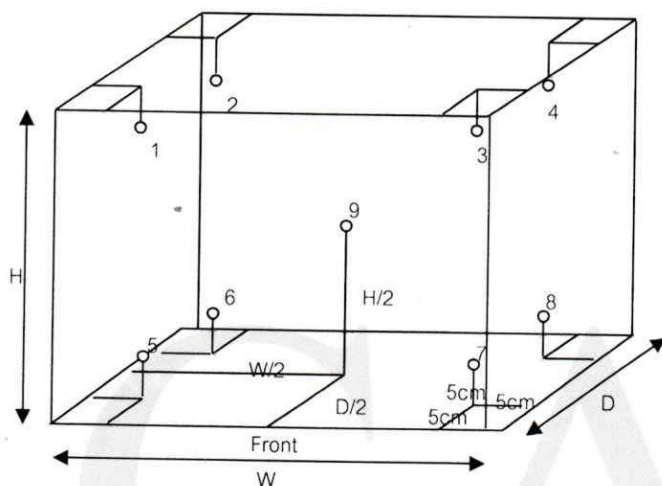

Inside of Chamber

W = 0.40 m

D = 0.33 m

H = 0.56 m

Capacity = 0.07 m³

Test Point (°C)	Setting Temperature (°C)	Indicating Temperature (°C)	Measured Temperature (°C) @ Sensor No.									Uncertainty (± °C)
			1	2	3	4	5	6	7	8	9	
85.0	85.0	85.0	85.3	85.1	85.2	85.1	85.0	84.8	84.7	84.7	85.0	0.64
104.0	104.0	104.0	104.3	104.2	104.3	104.1	104.0	103.8	103.7	103.7	104.0	0.69
180.0	180.0	180.0	179.5	179.9	180.0	179.8	179.9	179.4	178.9	179.4	180.0	0.95


Test Point (°C)	Setting Temperature (°C)	Indicating Temperature (°C)	Measured Uniformity (°C)	Measured Stability (°C)	Overall Variation (°C)
85.0	85.0	85.0	0.4	0.1	0.8
104.0	104.0	104.0	0.4	0.1	0.8
180.0	180.0	180.0	1.4	0.3	1.6

Remark The uncertainty is not combine uniformity of the air chamber

This result of calibration was found accurate as shown on date and place of calibration only.

This reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%

- o0o -





ID LINE : IEC17025



Certificate of Calibration

Certificate Number : SPR23110535-2

Page : 1 of 3

Customer : TNP ENVIRONMENT CO.,LTD.

332/173 Moo.3, Bang Rak Phatthana, Bang Bua Thong, Nonthaburi
11110

Equipment Name : Water Quality Meter

Manufacturer : Digicon

Model : WA-48SD

Serial Number : T.075714

ID. Number : TNP-LAB-46

Environmental Conditions

Ambient Temperature : $23^{\circ}\text{C} \pm 2^{\circ}\text{C}$

Relative Humidity : $50\% \pm 15\%$

Location of Calibration : In-Lab

Calibration Procedure : SP-CPC-04-01,
SP-CPC-04-02,

Method of Calibration : SP-CPC-04-11

Received Date : 30 Nov 2023

Calibration Date : 19 Jan 2024

Recommend Due Date : 19 Jan 2025

Date of Issue : 20 Jan 2024

This certifies that the above instrument was calibrated in compliance with the calibration system requirement of ISO/IEC 17025:2017 in accordance with reference procedure. Standards used to perform this calibration are certified by to NIST or equivalent, National metrology institute, Natural physical constants, consensus standards. The result reported herein apply only to the calibration of the item described above as received. Our decision rule is to contact the customer if the item pass and fail calibration when the results include the uncertainties and the customer must determine if the results meets their needs.

The calibration certificate shall not be reproduced except in full, without written approval of SP Metrology



ID LINE : IEC17025



Calibration Report

Certificate Number : SPR23110535-2

Page : 2 of 3

Reference Standards

Equipment Name	Model	Serial No.	Certificate No.	Due. Date
Zero Oxygen Solution	HI7040L	Lot S0027-23 _	21C31	21 Mar 2028
Standard pH Solution	PH016.L5	Lot No.882984	61267077	20 Mar 2024
Standard pH Solution	PH107.L5	Lot No.882985	61275614	13 Apr 2024
Standard pH Solution	PH020.L5	Lot No.882986	61268050	20 Mar 2024
Conductivity Standard 84 uS/cm	CS84M0S.L5	Lot No.882987	61247444	20 Mar 2024
Conductivity Standard 1413 uS/cm	CS1413M0S.L5	Lot No.882988	61267992	20 Mar 2024
Conductivity Standard 12880 uS/cm	CS1288P1S.L5	Lot No.882989	61247253	20 Mar 2024
Sodium Chloride Standard Solution	RM003461L25	Lot No.841770	97756699	23 Sep 2024

Traceability

This certification is traceable to the International System of Unit maintained at :

HANNA - Hanna Instruments (Thailand) Ltd.

C.P.A. Chem - ANAB#AT-1836 (ISO/IEC 17025:2017) and ANAB#AR-1835 (ISO/IEC 17034:2016)



ID LINE : IEC17025



Result of Calibration

Certificate No. : SPR23110535-2

Page : 3 of 3

pH Measurement @ 25 °C

Unit : pH

Standard Solution	UUC Reading	Error	Uncertainty (±)
4.008	3.98	-0.028	0.014
6.984	7.01	0.026	0.012
10.011	10.02	0.009	0.018

Conductivity Measurement @ 25 °C

Standard Solution	UUC Reading	Error	Uncertainty (±)
84 µS/cm	84.6 µS/cm	0.6 µS/cm	0.60 µS/cm
1.413 mS/cm	1.415 mS/cm	0.002 mS/cm	0.0082 mS/cm
12.88 mS/cm	12.89 mS/cm	0.01 mS/cm	0.075 mS/cm

* Dissolved Oxygen Permanance Test

Unit : mg/L

Actual Standard	UUC Reading	Error	Uncertainty (±)
0.0	0.0	0.0	0.13
5.0	4.9	-0.1	0.13
8.3	8.2	-0.1	0.13

Salinity Measurement

Unit : % Salinity

Actual Standard	UUC Reading	Error	Uncertainty (±)
5.0269	5.04	0.0131	0.015

Note:

The result of calibration was found accurate as show on date and place of calibration only.
This Certificate is not certified for any commercial transaction.
Calibration Marked (*) "Not ANAB Accredited " in this Certificate have been included for completeness.

Measurement Uncertainty

The reported uncertainty of measurement is the expanded uncertainty obtained by multiplying the standard uncertainty with the coverage factor $k = 2$, providing a level of confidence approximately 95%.

- End of Certificate -



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000-29 FAX. 0-2719-9484



Cert. No.: 23TM1099

Page : 1 of 3

Certificate of Calibration

Equipment :	Autoclave
Manufacturer :	Rexall
Model :	LS-2D
Serial No. :	04131
ID No. :	AUT-01
Submitted by :	Environment & Laboratory Co.,Ltd. 40 Soi Liangmueangnonthaburi 13, Talad Kwan, Mueang, Nonthaburi 11000
Location :	Room No. 205
Received Order :	12 July 2023
Calibration Date :	12 July 2023
Ambient Temperature :	(26 ± 10) °C
Relative Humidity :	(50 ± 30) %

Issue Date :

24 July 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

A 0056477



Equipment : Autoclave
Condition As-Received : Used Item
Reference : 2307-0094OC-7

Cert. No.: 23TM1099

Page : 2 of 3

Procedure Used :-

Calibration were conducted using in-house calibration procedure CP-OT03 according to direct measurement method with Data Acquisition which connected with Thermocouple Type T

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

<u>Instrument</u>	<u>Serial No.</u>	<u>Cert. No.</u>	<u>Traceable</u>	<u>Due Date</u>
1) Data Acquisition	MY41021843	22LM172	TPA	27 Dec 2023

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

4. This result of calibration covers laboratory autoclaves for the sterilization of goods and material which could be infected with organisms categorized as Hazard Group 1, 2 and 3**

(** = Categorization of pathogens according to hazard and categories of containment, second edition, 1990)


It does not cover autoclaves for use with material infect with organisms in Hazard Group 4, for which complete containment and sterilization of infected condensate is considered to be essential.

This result of calibration does not apply to sterilizers or disinfectors used for medical, dental, pharmaceutical or veterinary purposes which are directly concerned with patient care, or those used for fabrics subjected to sterilization which are required to be dry at the end of cycle.

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

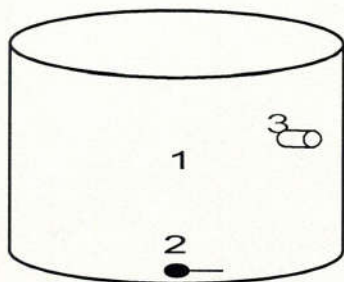

	Environmental		
	(°C)	(%R.H.)	(Volt)
Beginning of Calibration	30	55	220
Finished of Calibration	30	57	220

<u>Position</u>	<u>Description</u>	<u>Ref. Std. ID No.:</u>
1 =	Center of chamber	21-04TC-01
2 =	Temperature sensor	21-04TC-02
3 =	Exhaust port	21-04TC-03



Equipment : Autoclave
Condition As-Received : Used Item
Reference : 2307-0094OC-7
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source

Cert. No.: 23TM1099
Page : 3 of 3

Operating parameter Set : Temperature = 121 °C
Sterilization period = 15 minute

UUC* Setting (°C)	UUC* Reading (°C)	Position	Average* Standard Reading (°C)	Stability (± °C)	Pressure Reading (kg/cm ²)	Uncertainty (± °C)	Coverage Factor <i>k</i>
121	-	1	121.837	0.89	1.2	1.3	2
		2	121.869				
		3	121.875				

Average* : The average of 30 values in each position.

Stability : One-half of the greatest maximum difference of measured temperature at any one probe.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000-29 FAX. 0-2719-9484



Cert. No.: 23TM1100

Page : 1 of 3

Certificate of Calibration

Equipment : Incubator

Manufacturer : Memmert

Model : BM 500

Serial No. : D593.0342

ID No. : CHI-002

Submitted by : Environment & Laboratory Co.,Ltd.
40 Soi Liangmueangnonthaburi 13,
Talad Kwan, Mueang,
Nonthaburi 11000

Location : Room No. 204

Received Order : 12 July 2023

Calibration Date : 12 July 2023

Ambient Temperature : (26 ± 10) °C

Relative Humidity : (50 ± 30) %

Issue Date :

24 July 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

A 0056478



Equipment : Incubator
 Condition As-Received : Used Item
 Reference : 2307-0094OC-4

Cert. No.: 23TM1100
 Page : 2 of 3

Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD).

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Data Acquisition	MY41021843	22LM172	TPA	27 Dec 2023

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.


Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Not Available

Environment during calibration		
	Beginning	Finished
Temp. (°C)	22	22
REL.Humid. (%)	64	66
AC Supply (Volt)	222	221

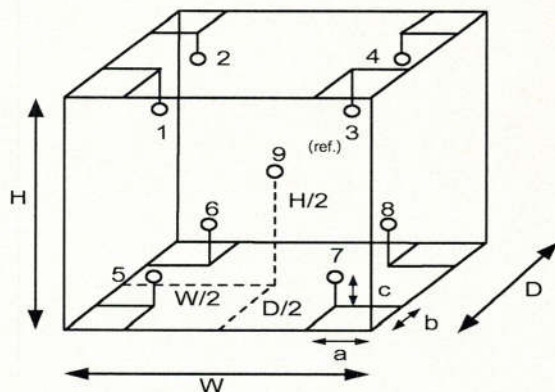

Position :	Ref. Std. ID No.:
1	18-04RTD-01
2	18-04RTD-02
3	18-04RTD-03
4	18-04RTD-04
5	18-04RTD-05
6	18-04RTD-06
7	18-04RTD-07
8	18-04RTD-08
9 (ref.)	18-04RTD-09

Probe Installation Details :

a = 5.0 cm
 b = 5.0 cm
 c = 5.0 cm

Dimension of Chamber :

D = 0.40 m
 W = 0.56 m
 H = 0.48 m
 Capacity = 0.11 m³



Equipment : Incubator
Condition As-Received : Used Item
Reference : 2307-0094OC-4
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Not Available

Cert. No.: 23TM1100

Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor <i>k</i>
35.0	35.0	35.0	0.040	0.27	0.48	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	
35.0	35.021	34.900	35.010	35.284	34.853	34.919	34.945	34.964	35.089	0.30

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000-29 FAX. 0-2719-9484



Cert. No.: 23TM1273

Page : 1 of 3

Certificate of Calibration

Equipment :	Incubator
Manufacturer :	Envilab-Intiligent
Model :	-
Serial No. :	-
ID No. :	CHI-005
Submitted by :	Environment & Laboratory Co.,Ltd. 40 Soi Liangmueangnonthaburi 13, Talad Kwan, Mueang, Nonthaburi 11000
Location :	Room No. 204
Received Order :	24 August 2023
Calibration Date :	24 August 2023
Ambient Temperature :	(26 ± 10) °C
Relative Humidity :	(50 ± 30) %

Issue Date :

29 August 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

A 0057741



Equipment : Incubator
Condition As-Received : Used Item
Reference : 2308-0600OC-1

Cert. No.: 23TM1273

Page : 2 of 3

Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD).

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Data Acquisition	MY44035217	22LM170	TPA	16 Dec 2023

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.


Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Not Available

Environment during calibration		
	Beginning	Finished
Temp. (°C)	23	23
REL.Humid. (%)	50	54
AC Supply (Volt)	220	220

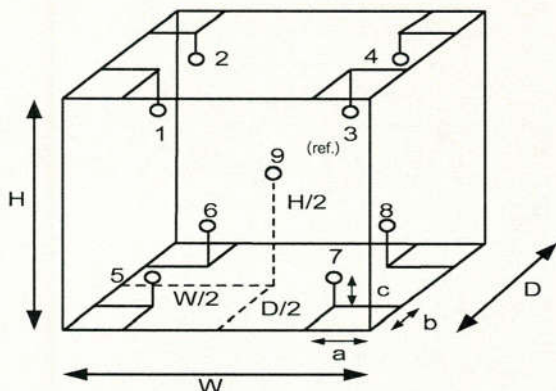

Position :	Ref. Std. ID No.:
1	18-10RTD-01
2	18-10RTD-02
3	18-10RTD-03
4	18-10RTD-04
5	18-10RTD-05
6	18-10RTD-06
7	18-10RTD-07
8	18-10RTD-08
9 (ref.)	18-10RTD-09

Probe Installation Details :

a = 5.0 cm
b = 5.0 cm
c = 5.0 cm

Dimension of Chamber :

D = 0.40 m
W = 0.70 m
H = 1.0 m
Capacity = 0.28 m³



Equipment : Incubator
Condition As-Received : Used Item
Reference : 2308-0600OC-1
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Not Available

Cert. No.: 23TM1273

Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor <i>k</i>
35.0	35.0	35.0	0.47	1.2	1.8	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	
35.0	34.805	34.737	34.701	34.435	34.724	34.783	35.228	35.604	34.816	0.71

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000-29 FAX. 0-2719-9484



Cert.No.: 23CHO420

Page.: 1 of 3

Certificate of Calibration

Equipment :	Spectrophotometer
Manufacturer :	Hach
Model :	DR 3900
Serial No. :	1988383
ID No. :	SPE-002
Condition As-Received:	Used Item
Received Date :	12 July 2023
Calibration Date :	12 July 2023
Reference :	2307-0094OC-11
Submitted by :	Environment & Laboratory Co.,Ltd. 40 Soi Liangmueangnonthaburi 13 Talad Kwan, Mueang, Nonthaburi 11000
Calibration Place :	Room No. 304
Ambient Temperature :	(26.7 - 26.9) °C (On-Site)
Relative Humidity :	(57.2 - 51.2) % (On-Site)
Calibration Procedure :	In - house method :

75-01

Issue Date :

24 July 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

A 0056486



Cert. No. : 23CHO420

Page : 2 of 3

Condition of calibration result

1. Reference Standard Material :

<u>Material</u>	<u>Serial No.</u>	<u>Certificate No.</u>	<u>Due date</u>
1. Absorbance Standard set	32587	100509	28 Mar 2024
2. Absorbance Standard set	32590	100508	28 Mar 2024
3. Absorbance Standard set	8331	105939	28 Sep 2024
4. Wavelength Standard set	29829	94776	02 Sep 2023
5. Wavelength Standard set	29829	94777	02 Sep 2023

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certificate is traceable to the International System of Unit maintained at :

- National Physical Laboratory (NPL), The United Kingdom of Great Britain and Northern Ireland
- National Institute of Standards and Technology (NIST), The United States of America

4. Spectral BandWidth : 5 nm

Scan Speed : - nm/min

Calibration Results : without adjustment

Wavelength Accuracy

Certified Values of Reference Material (nm)	UUC Reading (nm)	Uncertainty of Measurement (\pm nm)	Coverage Factor <i>k</i>
361.40	361	0.59	2.00
447.20	446	0.59	2.00
537.00	536	0.59	2.00
638.00	637	0.59	2.00
740.51	741	0.59	2.00
807.04	807	0.59	2.00



Cert. No. : 23CHO420

Page : 3 of 3

Calibration Results : without adjustment

Photometric Accuracy

Wavelength (nm)	Certified Values of Reference Material (Abs)	UUC Reading (Abs)	Uncertainty of Measurement (\pm Abs)	Coverage Factor <i>k</i>
350.0	Zero	0.000	0.0046	2.00
	0.4246	0.423	0.0061	2.00
	Zero	0.000	0.0050	2.00
	0.6385	0.633	0.0055	2.00
440.0	Zero	0.000	0.0028	2.00
	0.5607	0.560	0.0030	2.00
	0.7336	0.733	0.0030	2.00
	1.0636	1.063	0.0030	2.00
546.1	Zero	0.000	0.0028	2.00
	0.5224	0.522	0.0028	2.00
	0.6856	0.684	0.0029	2.00
	0.9937	0.992	0.0028	2.00
635.0	Zero	0.000	0.0028	2.00
	0.5397	0.538	0.0028	2.00
	0.6832	0.680	0.0028	2.00
	0.9886	0.985	0.0028	2.00

Remark

- Each individual filter is measured against the empty filter holder (blank) used to zero the spectrophotometer
- The Potassium Dichromate filled cells are measured against a Perchloric acid blank.

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k* , providing a level of confidence of approximately 95 %.

-o0o-





TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000-29 FAX. 0-2719-9484



Cert. No.: 23TM1173

Page : 1 of 3

Certificate of Calibration

Equipment :	Water Bath
Manufacturer :	Memmert
Model :	WB 22
Serial No. :	I505.0053
ID No. :	WAB-01
Submitted by :	Environment & Laboratory Co.,Ltd. 40 Soi Liangmueangnonthaburi 13, Talad Kwan, Mueang, Nonthaburi 11000
Location :	Room No. 303
Received Order :	12 July 2023
Calibration Date :	12 - 13 July 2023
Ambient Temperature :	(26 ± 10) °C
Relative Humidity :	(50 ± 30) %

Issue Date :

24 July 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

A 0056487



Equipment : Water Bath
Condition As-Received : Used Item
Reference : 2307-0094OC-3

Cert. No.: 23TM1173
Page : 2 of 3

Procedure Used :-

Calibration were conducted using in-house calibration procedure CP-OT04 according to direct measurement method with Data Acquisition which connected with Industrial Platinum Resistance Thermometer (IPRT).

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

<u>Instrument</u>	<u>Serial No.</u>	<u>Cert. No.</u>	<u>Traceable</u>	<u>Due Date</u>
1) Data Acquisition	MY44073381	23LM95	TPA	19 May 2024

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.


Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Heat transfer medium used : Water

	Environmental		AC Voltage Supply
	(°C)	(%R.H.)	
Beginning of Calibration	30	47	220
Finished of Calibration	31	50	221

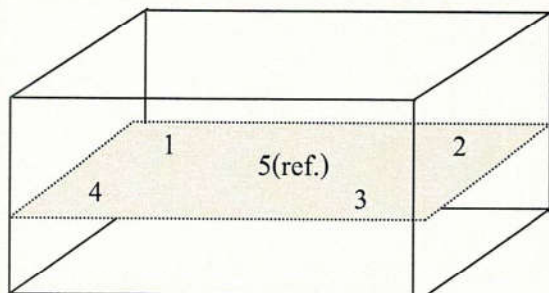

Front

Position :	Ref. Std. S/N.:
1	4803988-006
2	4803988-007
3	4804539-014
4	4804539-015
5(ref.)	4804539-016



Equipment : Water Bath
Condition As-Received : Used Item
Reference : 2307-0094OC-3
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source

Cert. No.: 23TM1173

Page : 3 of 3

Calibration point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Average* Standard Reading (°C)					Uncertainty (± °C)
			Position					
			1	2	3	4	5 (ref.)	
44.5	44.5	44.5	44.507	44.503	44.498	44.509	44.502	0.15
60.0	60.0	60.0	59.914	59.928	59.912	59.899	59.894	0.15

Calibration point (°C)	Uniformity (°C)	Stability (± °C)	Coverage Factor <i>k</i>
44.5	0.039	0.023	2
60.0	0.098	0.042	2

Average* : The average of 30 values in each position.

Uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Stability : One-half of the greatest maximum difference of measured temperature at any one probe.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity.

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-



CALIBRATION CERTIFICATE

Date of Issue Jun 28, 2024

Cert No. 24/2419

Site Calibration

Order No. 24060337

Customer Bureau Veritas AQ Lab (Thailand) Limited
111 Thailand Science Park
Moo 9 Paholyothin Rd., Klong 1, Klong Luang, Pathumthani 12120 Thailand

Place of Calibration Prepare Culture Media

Description Autoclave

Model VE-65

Serial No. 7202

ID.No. CHM000206

Date of Receipt Jun 26, 2024

Date of Calibration Jun 26, 2024

Environment

Temperature	(Min)	22.1	°C	(Max)	25.5	°C
Relative Humidity	(Min)	42.7	%RH	(Max)	56.3	%RH

Calibration Method

WI-22 : The reference thermometer was placing into the chamber and measurement based on BS-2646.
The temperature scale in use at this laboratory is the International Temperature Scale of 1990.



CALIBRATION CERTIFICATE


Date of Issue Jun 28, 2024

Cert No. 24/2419

Site Calibration

Order No. 24060337

Results (without adjustment)

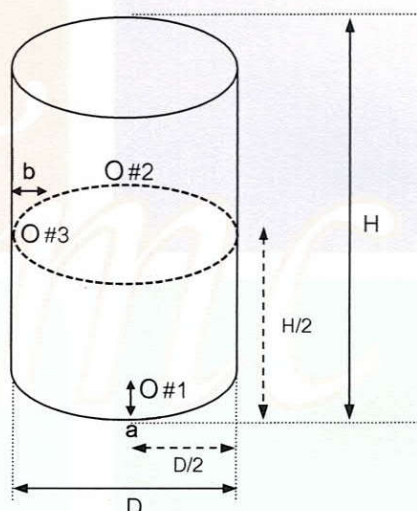

Position of reference thermometers were placed

Note.

- 1). Duration of the holding time was 45 minutes.
- 2). Stability - greatest one half of difference between max peak and min peak of each reference probe measured temperature obtained during the calibration interval.
- 3). Uniformity - the maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady state conditions. The reference sensor should preferably be located at the geometric center of the chamber.



CALIBRATION CERTIFICATE

Date of Issue Jun 28, 2024

Cert No. 24/2419

Site Calibration

Order No. 24060337

Results (without adjustment)

UUC Setting (°C)	UUC Reading		Reference Thermometer (°C)		Stability ±(°C)	Uniformity (°C)	Uncertainty ±(°C)
	(°C)	(KPa)					
115.0	115.3	179.6	Position 1	115.58	0.25	0.15	0.70
			Position 2	115.58			
			Position 3	115.48			

UUC Setting (°C)	UUC Reading		Reference Thermometer (°C)		Stability ±(°C)	Uniformity (°C)	Uncertainty ±(°C)
	(°C)	(KPa)					
121.0	121.4	211.0	Position 1	121.59	0.20	0.13	0.70
			Position 2	121.60			
			Position 3	121.52			


Standard

1) Standard Datalogger with RTD Model. HiTemp. 140-PT S/N. R32397, Certificate No. QR23-1962, Calibrated by Quality Reborn Co., Ltd., ONAC Calibration No. 0292. Due Date Aug 07, 2024.

2) Standard Datalogger with RTD Model. HiTemp. 140-PT S/N. R32398, Certificate No. QR23-1963, Calibrated by Quality Reborn Co., Ltd., ONAC Calibration No. 0292. Due Date Aug 07, 2024.

3) Standard Datalogger with RTD Model. HiTemp. 140-PT S/N. R32399, Certificate No. QR23-1964, Calibrated by Quality Reborn Co., Ltd., ONAC Calibration No. 0292. Due Date Aug 07, 2024.

This certificate is traceable to SI unit.

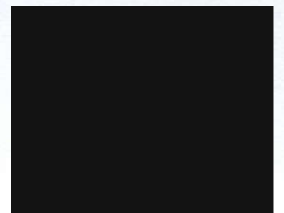



Thermology Co., Ltd.

96/177-96/178 Moo 6, T. La-harn, A. Bangbuathong, Nonthaburi 11110
Tel : 0 2191 6479 Fax : 0 2191 6480 website : www.thermology.co



CALIBRATION CERTIFICATE

Date of Issue Jun 28, 2024

Site Calibration

Cert No. 24/2419

Order No. 24060337

The stability and uniformity was taken into account in the measurement uncertainty stated.

The above results are valid exclusively for calibration samples as mentioned in the report.

The reported expanded uncertainty is based on a standard uncertainty multiplied by a coverage factor $k=2$, providing a level of confidence of approximately 95%. The uncertainty evaluation has been carried out in accordance with ONAC requirements.




Thermology

Autoclave

Cert.No. 24/2419

Model. VE-65 S/N. 7202 ID.No. CHM000206

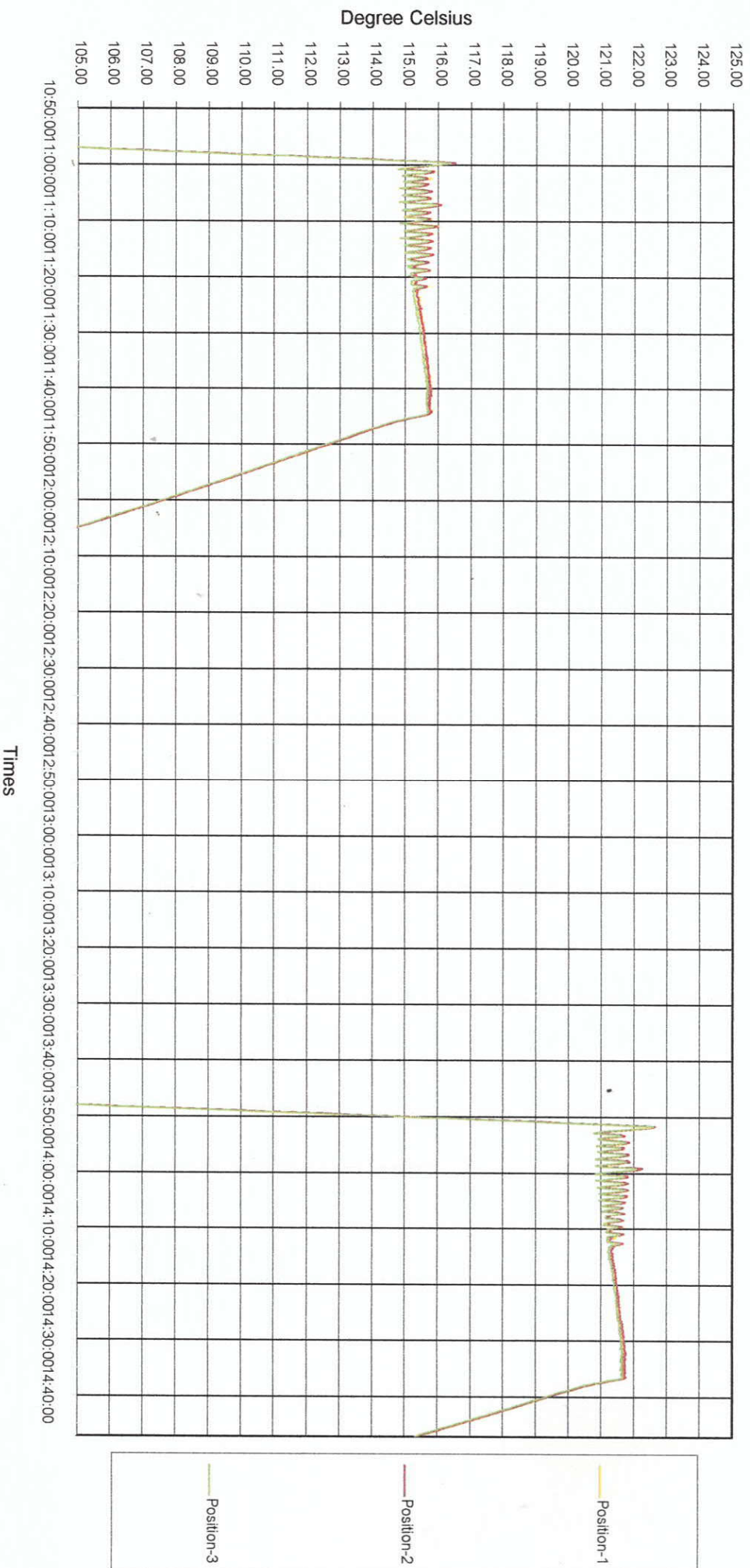



CALIBRATION CERTIFICATE

Date of Issue Jun 28, 2024

Cert No. 24/2418

Site Calibration

Order No. 24060337

Customer Bureau Veritas AQ Lab (Thailand) Limited
111 Thailand Science Park
Moo 9 Paholyothin Rd., Klong 1, Klong Luang, Pathumthani 12120 Thailand

Place of Calibration Incubation Room

Description Water Bath

Model SC100

Serial No. 0152187501160414

ID.No. CHM000205

Date of Receipt Jun 24, 2024

Date of Calibration Jun 24, 2024

Environment

Temperature	(Min)	22.8	°C	(Max)	25.2	°C
Relative Humidity	(Min)	44.1	%RH	(Max)	58.5	%RH
Line Voltage	(Min)	227.2	VAC	(Max)	229.6	VAC

Calibration Method

WI-18 : The reference thermometer was placed into the chamber and measurement was performed based on AS-2853.

The temperature scale in use at this laboratory is the International Temperature Scale of 1990.

Standard

1) Data Acquisition with Sensor Model 34972A S/N. MY49007789, Certificate No. QR24-0186, Calibrated by Quality Reborn Co., Ltd., ONAC Calibration No. 0202, Due Date Jan 23, 2025

This certificate is traceable to

This certificate is issued in accordance with the requirements of Thermology Laboratory and is valid only for the unit of measurement realised at corresponding national standard laboratory. This certificate may not be reproduced other than in full except with the prior written approval of laboratory.



CALIBRATION CERTIFICATE


Date of Issue Jun 28, 2024

Site Calibration

Cert No. 24/2418

Order No. 24060337

Results (without adjustment)

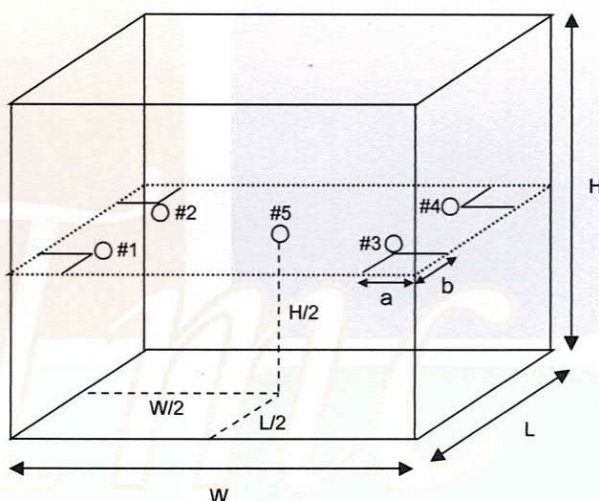

Position of reference thermometers were placed

Note.

- 1). Dimension (W x L x H) is 30 x 34 x 20 cm.
- 2). Stability - greatest one half of difference between max peak and min peak of each reference probe measured temperature obtained during the calibration interval.
- 3). Uniformity - the maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady state conditions. The reference sensor should preferably be located at the geometric center of the chamber.



CALIBRATION CERTIFICATE

Date of Issue Jun 28, 2024

Cert No. 24/2418

Site Calibration

Order No. 24060337


Results (without adjustment)

Cal Point (°C)	UUC Setting (°C)	UUC Reading (°C)	Reference Thermometer (°C)		Stability \pm (°C)	Uniformity (°C)	Uncertainty \pm (°C)
44.5	44.5	44.5	Position 1	44.490	0.020	0.043	0.13
			Position 2	44.494			
			Position 3	44.491			
			Position 4	44.499			
			Position 5	44.503			

The stability and uniformity was taken into account in the measurement uncertainty stated.

The above results are valid exclusively for calibration samples as mentioned in the report.


The reported expanded uncertainty is based on a standard uncertainty multiplied by a coverage factor $k=2$, providing a level of confidence of approximately 95%. The uncertainty evaluation has been carried out in accordance with ONAC requirements.

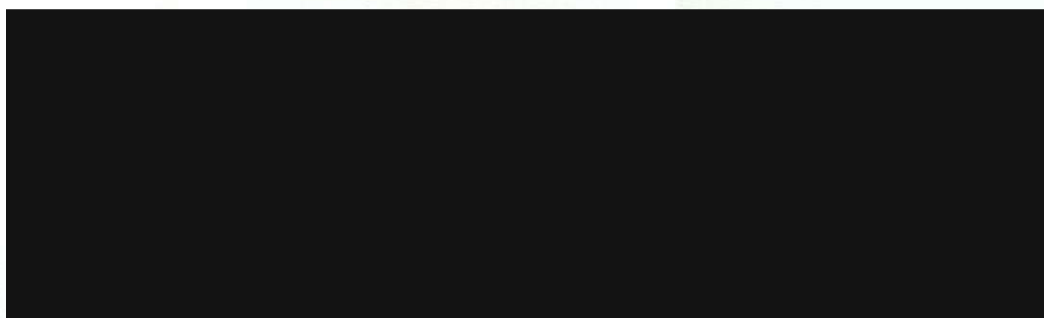

Water Bath

Cert.No. 24/2418

Model. SC100 S/N. 0152187501160414 ID.No. CHM000205

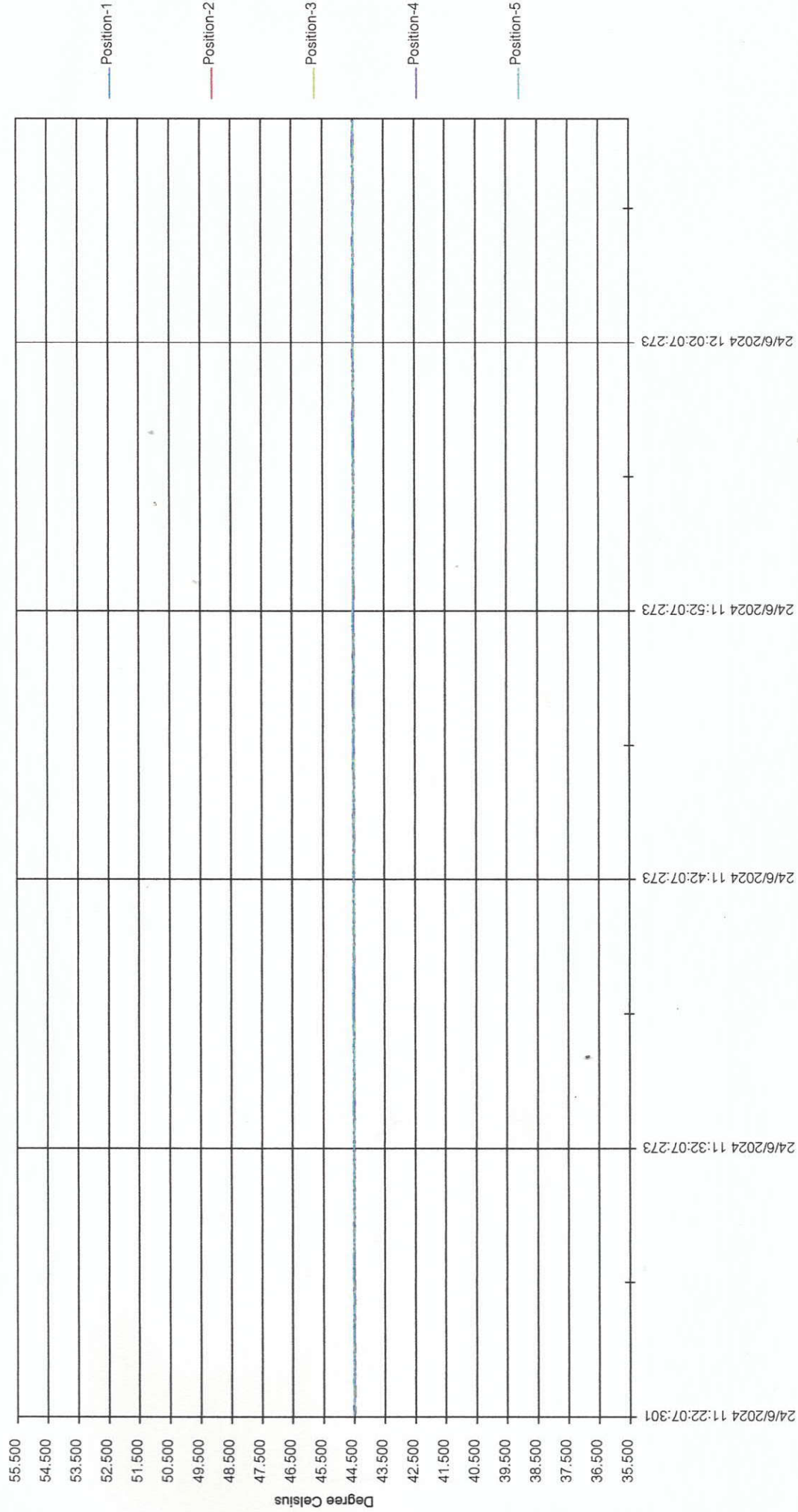



CALIBRATION CERTIFICATE

Date of Issue Jun 28, 2024

Cert No. 24/2415

Site Calibration

Order No. 24060337

Customer Bureau Veritas AQ Lab (Thailand) Limited
111 Thailand Science Park
Moo 9 Paholyothin Rd., Klong 1, Klong Luang, Pathumthani 12120 Thailand

Place of Calibration Incubation Room

Description Incubator

Model IN110

Serial No. D415.0797

ID.No. CHM000181

Date of Receipt Jun 24, 2024

Date of Calibration Jun 24, 2024

Environment

Temperature (Min) 22.8 °C (Max) 25.2 °C

Relative Humidity (Min) 44.1 %RH (Max) 58.5 %RH

Calibration Method

WI-17 : The reference thermometer was placed into the chamber and measurement was performed based on AS-2853.

The temperature scale in use at this laboratory is the International Temperature Scale of 1990.

Standard

1) Data Acquisition with Sensor Model 34972A S/N. MY49010059, Certificate No. QR24-0874, Calibrated by Quality Reborn Co., Ltd., ONAC Calibration No. 0292. Due Date Apr 24, 2025.

This certificate is traceable to SI unit.



CALIBRATION CERTIFICATE


Date of Issue Jun 28, 2024

Site Calibration

Cert No. 24/2415

Order No. 24060337

Results (without adjustment)

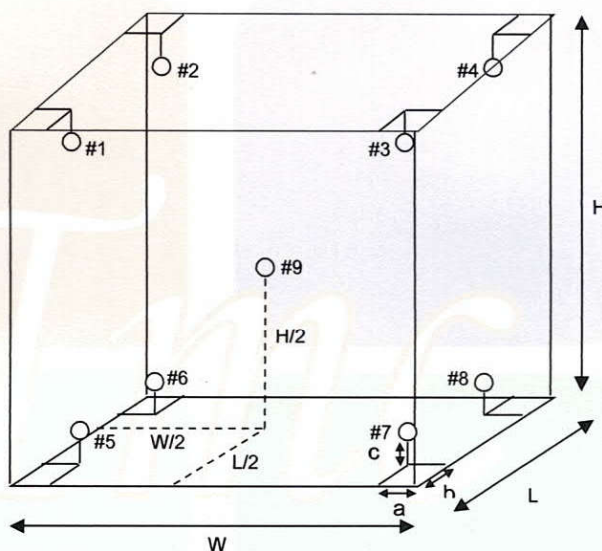

Position of reference thermometers were placed

Note.

- 1). Dimension ($W \times L \times H$) is 56 x 40 x 48 cm.
- 2). Stability - greatest one half of difference between max peak and min peak of each reference probe measured temperature obtained during the calibration interval.
- 3). Uniformity - the maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady state conditions. The reference sensor should preferably be located at the geometric center of the chamber.



CALIBRATION CERTIFICATE

Date of Issue Jun 28, 2024

Cert No. 24/2415

Site Calibration

Order No. 24060337

Results (without adjustment)

Cal Point (°C)	UUC Setting (°C)	UUC Reading (°C)	Reference Thermometer (°C)		Stability \pm (°C)	Uniformity (°C)	Uncertainty \pm (°C)
35.0	35.0	35.0	Position 1	35.138	0.067	0.253	0.30
			Position 2	35.099			
			Position 3	35.075			
			Position 4	35.187			
			Position 5	35.173			
			Position 6	34.988			
			Position 7	34.878			
			Position 8	34.965			
			Position 9	34.970			


The stability and uniformity was taken into account in the measurement uncertainty stated.

The above results are valid exclusively for calibration samples as mentioned in the report.

The reported expanded uncertainty is based on a standard uncertainty multiplied by a coverage factor $k=2$, providing a level of confidence of approximately 95%. The uncertainty evaluation has been carried out in accordance with ONAC requirements.

Incubator
Model: IN110 S/N: D415.0797 ID.No: CHM000181

Cert.No. 24/2415

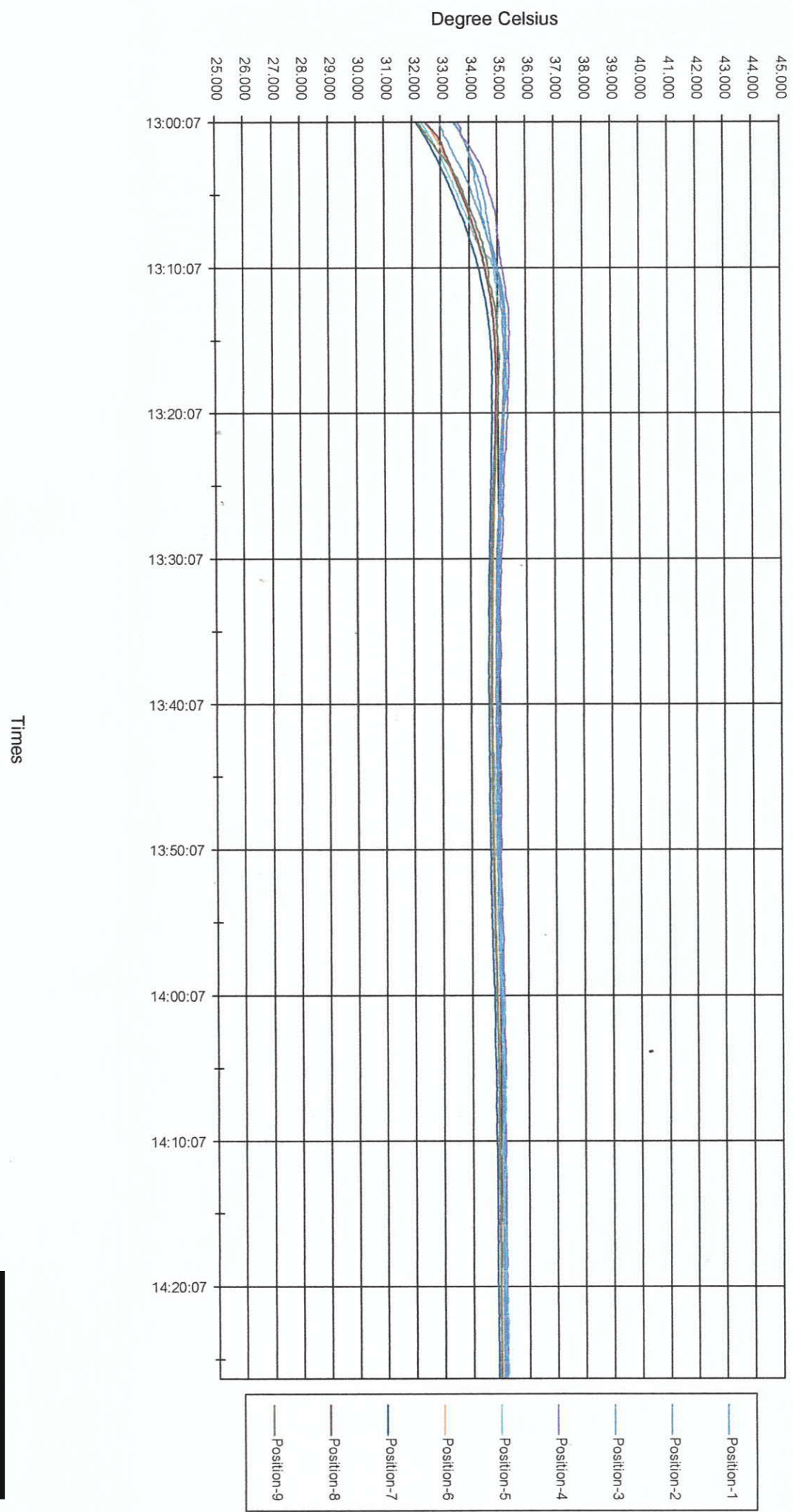



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000-29 FAX. 0-2719-9484



Cert.No.: 23CHO565

Page.: 1 of 3

Certificate of Calibration

Equipment : Spectrophotometer
Manufacturer : Merck
Model : Prove 100
Serial No. : 1904113676
ID No. : S2019025
Condition As-Received: Used Item
Received Date : 21 September 2023
Calibration Date : 21 September 2023
Reference : 2309-0483OC-2
Submitted by : SGS (Thailand) Limited
1/209, 1/211 Moo 1, Ban Chang,
Ban Chang, Rayong 21130

Calibration Place : Spectrophotometry Lab
Ambient Temperature : (23.8 - 20.9) °C (On-Site)
Relative Humidity : (50.1 - 50.2) % (On-Site)
Calibration Procedure : In - house method :
GB 6684 based on ASTM E 275-01

Issue Date : 20 September 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

A 0058836



Cert. No. : 23CHO565

Page : 2 of 3

Condition of calibration result

1. Reference Standard Material :

<u>Material</u>	<u>Serial No.</u>	<u>Certificate No.</u>	<u>Due date</u>
1. Absorbance Standard set	39130	106269	10 Oct 2024
2. Wavelength Standard set	36730	98330	19 Jan 2024
3. Wavelength Standard set	36730	98331	19 Jan 2024

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certificate is traceable to the International System of Unit maintained through :
- Starna Scientific Ltd.

4. Spectral BandWidth : 4 nm
Scan Speed : - nm/min

Calibration Results : without adjustment

Wavelength Accuracy

Certified Values of Reference Material (nm)	UUC Reading (nm)	Uncertainty of Measurement (\pm nm)	Coverage Factor <i>k</i>
418.48	418.5	0.15	2.00
513.70	513.1	0.14	2.00
536.90	536.3	0.14	2.00
637.94	637.6	0.14	2.00
879.70	878.8	0.15	2.00



Cert. No. : 23CHO565

Page : 3 of 3

Calibration Results : without adjustment

Photometric Accuracy

Wavelength (nm)	Certified Values of Reference Material (Abs)	UUC Reading (Abs)	Uncertainty of Measurement (\pm Abs)	Coverage Factor <i>k</i>
440.0	Zero	0.000	0.0028	2.00
	0.5645	0.563	0.0028	2.00
	0.6988	0.698	0.0028	2.00
	1.0017	1.001	0.0028	2.00
546.1	Zero	0.000	0.0028	2.00
	0.5281	0.527	0.0028	2.00
	0.6962	0.696	0.0028	2.00
	0.9984	0.998	0.0028	2.00
635.0	Zero	0.000	0.0028	2.00
	0.5699	0.569	0.0028	2.00
	0.7606	0.760	0.0028	2.00
	1.0927	1.092	0.0028	2.00

Remark

- Each individual filter is measured against the empty filter holder (blank) used to zero the spectrophotometer

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k* , providing a level of confidence of approximately 95 %.

-o0o-

Verification COD Reactor

Equipment Name	Dri-Block Heater-Digital	Temperature Ver	150±2 °C
Serial No.	000827-A	Model	DB 200/3
Reference Standard	Thermocouple Type K	Certificate No.	21/4272
Calibration Date	10/03/2023	Next Cal. Date	10/03/2024

Left											
Hole 1				Hole 2				Hole 3			
NO.	Result			NO.	Result			NO.	Result		
	temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.
1	151.1	-0.36	150.7	1	150.6	-0.36	150.2	1	151.4	-0.36	151.0
2	150.8	-0.36	150.4	2	151.7	-0.36	151.3	2	151.3	-0.36	150.9
3	151.2	-0.36	150.8	3	151.1	-0.36	150.7	3	151.7	-0.36	151.3
		Mean	150.67			Mean	150.77			Mean	151.11
		SD	0.208			SD	0.551			SD	0.208
		%RSD	0.138			%RSD	0.365			%RSD	0.138
Hole 4				Hole 5				Hole 6			
NO.	Result			NO.	Result			NO.	Result		
	temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.
1	151.7	-0.36	151.3	1	150.5	-0.36	150.1	1	151.5	-0.36	151.1
2	151.6	-0.36	151.2	2	151.3	-0.36	150.9	2	151.4	-0.36	151.0
3	151.5	-0.36	151.1	3	150.6	-0.36	150.2	3	150.5	-0.36	150.1
		Mean	151.24			Mean	150.44			Mean	150.77
		SD	0.100			SD	0.436			SD	0.551
		%RSD	0.066			%RSD	0.290			%RSD	0.365
Hole 7				Hole 8				Hole 9			
NO.	Result			NO.	Result			NO.	Result		
	temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.
1	151.3	-0.36	150.9	1	151.7	-0.36	151.3	1	150.5	-0.36	150.1
2	151.0	-0.36	150.6	2	150.5	-0.36	150.1	2	151.2	-0.36	150.8
3	151.3	-0.36	150.9	3	151.4	-0.36	151.0	3	150.8	-0.36	150.4
		Mean	150.84			Mean	150.84			Mean	150.47
		SD	0.173			SD	0.624			SD	0.351
		%RSD	0.115			%RSD	0.414			%RSD	0.233
Hole 10				Hole 11				Hole 12			
NO.	Result			NO.	Result			NO.	Result		
	temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.
1	151.6	-0.36	151.2	1	151.5	-0.36	151.1	1	150.7	-0.36	150.3
2	151.6	-0.36	151.2	2	151.2	-0.36	150.8	2	151.6	-0.36	151.2
3	150.8	-0.36	150.4	3	151.5	-0.36	151.1	3	151.2	-0.36	150.8
		Mean	150.97			Mean	151.04			Mean	150.81
		SD	0.462			SD	0.173			SD	0.451
		%RSD	0.306			%RSD	0.115			%RSD	0.299

Verification COD Reactor

Equipment Name	Dri-Block Heater Digital	Temperature Ver	150±2 °C
Serial No.	000827-A	Model	DB 200/3
Reference Standard	Thermocouple Type K	Certificate No.	21/4272
Calibration Date	10/03/2023	Next Cal. Date	10/03/2024

Middle

Hole 1				Hole 2				Hole 3			
NO.	Result			NO.	Result			NO.	Result		
	temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.
1	151.2	-0.36	150.8	1	150.7	-0.36	150.3	1	151.1	-0.36	150.7
2	151.5	-0.36	151.1	2	151.7	-0.36	151.3	2	151.6	-0.36	151.2
3	151.6	-0.36	151.2	3	150.8	-0.36	150.4	3	150.9	-0.36	150.5
		Mean	151.07			Mean	150.71			Mean	150.84
		SD	0.208			SD	0.551			SD	0.361
		%RSD	0.138			%RSD	0.365			%RSD	0.239

Hole 4				Hole 5				Hole 6			
NO.	Result			NO.	Result			NO.	Result		
	temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.
1	151.5	-0.36	151.1	1	151.6	-0.36	151.2	1	150.5	-0.36	150.1
2	150.7	-0.36	150.3	2	151.2	-0.36	150.8	2	150.6	-0.36	150.2
3	151.6	-0.36	151.2	3	151.2	-0.36	150.8	3	151.2	-0.36	150.8
		Mean	150.91			Mean	150.97			Mean	150.41
		SD	0.493			SD	0.231			SD	0.379
		%RSD	0.327			%RSD	0.153			%RSD	0.252

Hole 7				Hole 8				Hole 9			
NO.	Result			NO.	Result			NO.	Result		
	temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.
1	151.5	-0.36	151.1	1	151.5	-0.36	151.1	1	151.2	-0.36	150.8
2	150.5	-0.36	150.1	2	151.7	-0.36	151.3	2	150.9	-0.36	150.5
3	150.8	-0.36	150.4	3	151.1	-0.36	150.7	3	151.4	-0.36	151.0
		Mean	150.57			Mean	151.07			Mean	150.81
		SD	0.513			SD	0.306			SD	0.252
		%RSD	0.341			%RSD	0.202			%RSD	0.167

Hole 10				Hole 11				Hole 12			
NO.	Result			NO.	Result			NO.	Result		
	temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.
1	151	-0.36	150.6	1	151.5	-0.36	151.1	1	151.2	-0.36	150.8
2	150.6	-0.36	150.2	2	150.6	-0.36	150.2	2	150.5	-0.36	150.1
3	151.5	-0.36	151.1	3	151.0	-0.36	150.6	3	151.1	-0.36	150.7
		Mean	150.67			Mean	150.67			Mean	150.57
		SD	0.451			SD	0.451			SD	0.379
		%RSD	0.299			%RSD	0.299			%RSD	0.251

Verification COD Reactor

Equipment Name Dri-Block Heater Digital
 Serial No. 000827/A
 Reference Standard Thermocouple Type K
 Calibration Date 10/03/2023

Temperature Ver 150±2 °C
 Model DB 200/3
 Certificate No. 21/4272
 Next Cal. Date 10/03/2024

Right

Hole 1				Hole 2				Hole 3			
NO.	Result			NO.	Result			NO.	Result		
	temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.
1	149.9	-0.36	149.5	1	150.3	-0.36	149.9	1	151.0	-0.36	150.6
2	151.1	-0.36	150.7	2	151.0	-0.36	150.6	2	151.0	-0.36	150.6
3	150.9	-0.36	150.5	3	149.9	-0.36	149.5	3	150.4	-0.36	150.0
		Mean	150.27			Mean	150.04			Mean	150.44
		SD	0.643			SD	0.557			SD	0.346
		%RSD	0.428			%RSD	0.371			%RSD	0.230

Hole 4				Hole 5				Hole 6			
NO.	Result			NO.	Result			NO.	Result		
	temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.
1	150.8	-0.36	150.4	1	150.0	-0.36	149.6	1	150.5	-0.36	150.1
2	151.0	-0.36	150.6	2	150.0	-0.36	149.6	2	150.8	-0.36	150.4
3	150.9	-0.36	150.5	3	150.7	-0.36	150.3	3	149.8	-0.36	149.4
		Mean	150.54			Mean	149.87			Mean	150.01
		SD	0.100			SD	0.404			SD	0.513
		%RSD	0.066			%RSD	0.270			%RSD	0.342

Hole 7				Hole 8				Hole 9			
NO.	Result			NO.	Result			NO.	Result		
	temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.
1	150.8	-0.36	150.4	1	151.1	-0.36	150.7	1	150.2	-0.36	149.8
2	150.9	-0.36	150.5	2	150.7	-0.36	150.3	2	150.2	-0.36	149.8
3	151.0	-0.36	150.6	3	151.1	-0.36	150.7	3	149.9	-0.36	149.5
		Mean	150.54			Mean	150.61			Mean	149.74
		SD	0.100			SD	0.231			SD	0.173
		%RSD	0.066			%RSD	0.153			%RSD	0.116

Hole 10				Hole 11				Hole 12			
NO.	Result			NO.	Result			NO.	Result		
	temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.		temp. °C	Corr.	temp+Corr.
1	150.6	-0.36	150.2	1	150.5	-0.36	150.1	1	150.9	-0.36	150.5
2	150.5	-0.36	150.1	2	150.9	-0.36	150.5	2	150.0	-0.36	149.6
3	149.9	-0.36	149.5	3	151.1	-0.36	150.7	3	150.5	-0.36	150.1
		Mean	149.97			Mean	150.47			Mean	150.11
		SD	0.379			SD	0.306			SD	0.451
		%RSD	0.252			%RSD	0.203			%RSD	0.300

สรุปผลการ Verify

ปรับอุณหภูมิ 154.0 °C แต่ค่าควบคุมมีอุณหภูมิ 150 ± 2 °C ทุกครั้ง



ARCHEMICA

Certificate of Calibration

Aquion RFIC : Anion System (ID#1054)

**This certificate is to verify that instrument below are calibrated
by**

Archemica Lab Co.,Ltd.

Aquion


S/N 220380025

AS-DV

S/N 2203880170

For

SGS (Thailand) Limited (Rayong Branch).

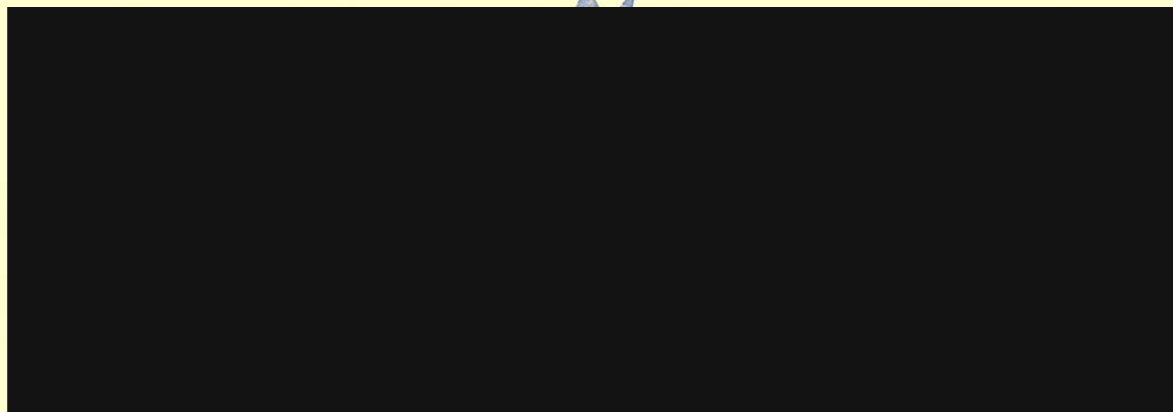

**MAINTENANCE AND TEST CERTIFICATE MODEL****Avio220 Max**

Date Tested: January 10, 2024
Recommendation Recertification
Period 6 Months
Recertification Due: July 10, 2024
Date Last Certified: N/A
Visit Number: 10E2 W
PerkinElmer Phone:
PerkinElmer Fax:

CONFIGURATION TESTED

MODEL	SERIAL NUMBER	SOFTWARE
AVIO200 MAX	M79S2304111	Syngristix V5
TESTED EQUIPMENT	CALIBRATION NUMBER	EXPIRATION
IPV Methods		
TEST STANDARD USED	PART NUMBER	EXPIRATION DATE
Multielement Standard	N069-1579	30-Dec-24
Instrument Cal. STD4	N930-0221	30-Nov-24
CUSTOMER SUPPLIED	COMMENTS	CUSTOMER INITIALS
2 % HNO3		
10 % HNO3		



MAINTENANCE AND TEST CERTIFICATE MODEL
Avio220 Max

SERIAL NUMBER M79S2304111**DATE TESTED** January 10, 2024**1. MECHANICAL CHECKS**

A. Inspect and clean all fans and filters.

☐ OK

B. Inspect and replace as necessary, all torch components including the RF coil.

☐ OK

C. Inspect all tubing for sign of clacking or leaking.

☐ OK

D. Adjust water and gas pressure regulator settings.

☐ OK

E. Inspect and leak check pneumatics drawers.

☐ OK

F. Clean the exterior of the instrument.

☐ OK**2. OPTICAL CHECKS**

A. Inspect and clean all optical components.

☐ OK

B. As required, check and replace all purgebfilters.

☐ OK

C. Recheck optical alignment.

☐ OK**3. COOLING SYSTEM CHECKS**

A. Perform preventive maintenance on chiller.

☐ OK

B. Flush out the chiller every year.

☐ OK**4. PERFORMANCE CHECKS**

A. Torch View Alignment.

☐ OK

B. Wavelength Calibration.

☐ OK



MAINTENANCE AND TEST CERTIFICATE MODEL

Avio220 Max

SERIAL NUMBER <u>M79S2304111</u>		DATE TESTED	<u>January 10, 2024</u>
PARAMETER	SPECIFICATION		FINAL VALUE
Spectral Resolution : UV			
As 193.696 nm	≤ 0.009 nm		<u>0.00864</u> nm
Ni 231.604 nm	≤ 0.011 nm		<u>0.01009</u> nm
Ni 341.476 nm	≤ 0.015 nm		<u>0.01169</u> nm
Spectral Resolution : VIS			
Ba 455.403 nm	≤ 0.020 nm		<u>0.01776</u> nm
Precision			
Zn 206.200 nm	% RSD ≤ 1.0 %		<u>0.28</u> %
Mg 280.271 nm	% RSD ≤ 1.0 %		<u>0.73</u> %
Mg 285.213 nm	% RSD ≤ 1.0 %		<u>0.61</u> %
Ba 455.403 nm	% RSD ≤ 1.0 %		<u>0.54</u> %
Detection Limits : Axial			
Tl 190.801 nm	3(sd)		<u>1.52</u> ppb
As 193.696 nm	3(sd)		<u>1.4</u> ppb
Se 196.026 nm	3(sd)		<u>1.53</u> ppb
Pb 220.353 nm	3(sd)		<u>1.72</u> ppb
Detection Limits : Radial			
As 193.696 nm	3(sd)		<u>1.69</u> ppb
Zn 213.857 nm	3(sd)		<u>0.42</u> ppb
Mn 257.610 nm	3(sd)		<u>0.1</u> ppb
La 379.478 nm	3(sd)		<u>0.61</u> ppb
Ba 455.403 nm	3(sd)		<u>0.13</u> ppb
Ba 493.408 nm	3(sd)		<u>0.1</u> ppb
BEC : Axial (IB X 1000)/(IS-IB)			
Mn 257.610 nm	≤ 30 ppb		<u>7.83</u> ppb
BEC : Radial (IB X 1000)/(IS-IB)			
Mn 257.610 nm	≤ 30 ppb		<u>17.57</u> ppb



MAINTENANCE AND TEST CERTIFICATE MODEL
Avio220 Max

SERIAL NUMBER M79S2304111

DATE TESTED January 10, 2024

Remarks :

Commissioning follow as commissioning performance sheets.

This is to certify that the above tests have been performed and the configuration tested

☒

meets

☐

does not meet

the PerkinElmer Specifications listed on this certificate.

ที่ อก ๐๓๑๐(๑)



กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๒ ๓ กุมภาพันธ์ ๒๕๖๗

เรื่อง ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

เรียน กรรมการผู้จัดการ บริษัท ทีเอ็นพี เอ็นไวรอนเมนท์ จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๑๙ ธันวาคม ๒๕๖๖

สิ่งที่ส่งมาด้วย เอกสารแนบท้ายหนังสือต่ออายุรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
บริษัท ทีเอ็นพี เอ็นไวรอนเมนท์ จำกัด จำนวน ๒ แผ่น

ตามหนังสือที่อ้างถึง บริษัท ทีเอ็นพี เอ็นไวรอนเมนท์ จำกัด ขอต่ออายุหนังสือรับขึ้นทะเบียน
ห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ว-๓๑๘ สถานที่ตั้งเลขที่ ๓๓๒/๑๗๓ หมู่ที่ ๓ ตำบลบางรักพัฒนา
อำเภอบางบัวทอง จังหวัดนนทบุรี ต่อกรมโรงงานอุตสาหกรรม นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้บริษัท ทีเอ็นพี เอ็นไวรอนเมนท์ จำกัด ต่ออายุหนังสือ
รับขึ้น



ค. ขอบข่ายสารมลพิษที่ได้รับขึ้นทะเบียนให้วิเคราะห์ในน้ำเสีย น้ำใต้ดิน อากาศเสีย และ
สิ่งปฏิกูลหรือวัสดุที่ไม่ใช้แล้ว ตามสิ่งที่ส่งมาด้วย

หนังสือฉบับนี้จะหมดอายุในวันที่ ๑๔ มกราคม ๒๕๗๐ หากประสงค์จะต่ออายุหนังสือ
รับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน สามารถยื่นคำขอผ่านระบบอิเล็กทรอนิกส์ได้ที่หน้าเว็บไซต์
กรมโรงงานอุตสาหกรรม

จึงเรียนมาเพื่อทราบ



กองวิจัยและเตือนภัยมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕

โทรสาร ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๙๙

ไปรษณีย์อิเล็กทรอนิกส์ saraban@diw.mail.go.th



เอกสารแนบท้ายหนังสือต่ออายุรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท ทีเอ็นพี เอ็นไวรอนเม้นท์ จำกัด

เลขทะเบียน ว-๓๑๘

ที่ อก ๐๓๑๐(๑)/

ลงวันที่ ๒๓ กุมภาพันธ์ ๒๕๖๗

ขอขยายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๑๑ รายการ

น้ำเสีย จำนวน 7 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Biochemical Oxygen Demand	1) 5-Day BOD Test, Azide Modification Method ^[2] 2) 5-Day BOD Test, Membrane Electrode Method ^[2]
2	Oil & Grease	Liquid-Liquid, Partition-Gravimetric Method ^[2]
3	pH	Electrometric Method ^[2]
4	Sulfide	Iodometric Method ^[2]
5	Temperature	Laboratory and Field Methods ^[2]
6	Total Dissolved Solids	Dried at 180 °C ^[2]
7	Total Suspended Solids	Dried from 103 to 105 °C ^[2]

น้ำใต้ดิน จำนวน 1 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	pH	Electrometric Method ^[2]

อากาศเสีย (ปล่องระบาย) จำนวน 2 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Opacity	Ringelmann's Method ^[1]
2	Total Suspended Particulate	Isokinetic Sampling, Gravimetric Method ^[3]

สิ่งปฏิกูลหรือวัสดุที่ไม่ใช้แล้ว จำนวน 1 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	pH	Electrometric Method ^[4,5]

เอกสารอ้างอิง...

เอกสารอ้างอิง

1. กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2549. เรื่อง กำหนดค่าปริมาณเขม่าควันที่เจือปนในอากาศที่ระบายออกจากปล่องของหม้อน้ำโรงสีข้าวที่ใช้แก๊สเป็นเชื้อเพลิง.

ราชกิจจานุเบกษา. 4 ธันวาคม 2549. เล่มที่ 123 ตอนพิเศษ 125ง.

2. APHA, AWWA, WEF. **Standard Methods for the Examination of Water and Wastewater.** 24th ed. Washington, DC: APHA, 2023.

3. United States Environmental Protection Agency. **Standards of Performance for New Stationary Sources.** 40 CFR 60. Appendix A, 2023.

4. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. **pH Electrometric Measurement. SW-846 Method 9040C,** 2004.

5. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. **Soil and Waste pH. SW-846 Method 9045D,** 2004.

ที่ อก ๐๓๑๐(๑)/



กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๒๑ มิถุนายน ๒๕๖๔

เรื่อง ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

เรียน กรรมการผู้จัดการ บริษัท เอ็นไวรอนเม้นท์ แอนด์ แลบลอราตอรี จำกัด

อ้างถึง ๑. คำขอขึ้นทะเบียน/ต่ออายุหนังสืออนุญาตขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

ลงวันที่ ๙ เมษายน ๒๕๖๔

๒. หนังสือบริษัท เอ็นไวรอนเม้นท์ แอนด์ แลบลอราตอรี จำกัด เลขที่ EN ๐๐๓/๒๕๖๔

ลงวันที่ ๙ เมษายน ๒๕๖๔

๓. หนังสือบริษัท เอ็นไวรอนเม้นท์ แอนด์ แลบลอราตอรี จำกัด เลขที่ EN ๐๐๕/๒๕๖๔

ลงวันที่ ๒๘ เมษายน ๒๕๖๔

สิ่งที่ส่งมาด้วย เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท เอ็นไวรอนเม้นท์ แอนด์ แลบลอราตอรี จำกัด จำนวน ๒ แผ่น

ตามหนังสือที่อ้างถึง ๑ ๒ และ ๓ บริษัท เอ็นไวรอนเม้นท์ แอนด์ แลบลอราตอรี จำกัด ขอต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ว-๐๒๙ สถานที่ตั้งเลขที่ ๔๐ ซอยเลี้ยวเมืองนนทบุรี ๑๓ ตำบลตลาดขวัญ อำเภอเมืองนนทบุรี จังหวัดนนทบุรี ต่อกรมโรงงานอุตสาหกรรม นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้บริษัท เอ็นไวรอนเม้นท์ แอนด์ แลบลอราตอรี จำกัด ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน โดยมีองค์ประกอบดังนี้


ตามสิ่งที่ส่งมา

หนังสือ...

หนังสือฉบับนี้จะหมดอายุในวันที่ ๑๘ เมษายน ๒๕๖๗ หากประสงค์จะต่ออายุหนังสือ
รับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ให้ยื่นคำขอต่ออายุพร้อมเอกสารประกอบคำขอต่อ
กรมโรงงานอุตสาหกรรมภายใน ๓๐ วัน ก่อนวันสิ้นอายุของหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
ซึ่งคำขอต่ออายุดังกล่าวขอรับได้ที่กรมโรงงานอุตสาหกรรม

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

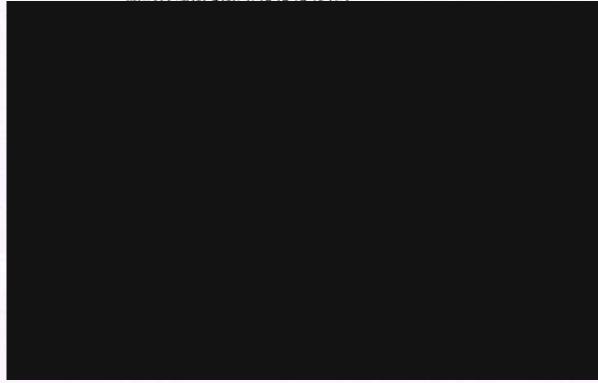

กองวิจัยและเตือนภัยมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๒๐๒ ๔๐๐๒ ๐ ๒๒๐๒ ๔๑๔๖

โทรสาร ๐ ๒๓๕๔ ๓๔๑๕

เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท เอ็นไวรอนเม้นท์ แอนด์ แลบลอราตอรี จำกัด

เลขทะเบียน ว-๐๒๙

ที่ อก ๐๓๑๐(๑)/

ลงวันที่ ๒๑ มิถุนายน ๒๕๖๕

ขอขยายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๒๕ รายการ

น้ำเสีย จำนวน 25 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Arsenic	Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[2]
2	Biochemical Oxygen Demand	5-Day BOD Test, Azide Modification Method ^[2]
3	Cadmium	Digestion, Direct Air-Acetylene Flame Method ^[2]
4	Chemical Oxygen Demand	Closed Reflux, Titrimetric Method ^[2]
5	Color	ADMI Weighted-Ordinate Spectrophotometric Method ^[2]
6	Copper	Digestion, Direct Air-Acetylene Flame Method ^[2]
7	Cyanide	Distillation, Colorimetric Method ^[2]
8	Formaldehyde	Distillation, Colorimetric Method ^[1]
9	Free Chlorine	Iodometric Method ^[2]
10	Hexavalent Chromium	Colorimetric Method ^[2]
11	Lead	Digestion, Direct Air-Acetylene Flame Method ^[2]
12	Manganese	Digestion, Direct Air-Acetylene Flame Method ^[2]
13	Mercury	Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[2]
14	Nickel	Digestion, Direct Air-Acetylene Flame Method ^[2]
15	Oil & Grease	1) Liquid-Liquid, Partition-Gravimetric Method ^[2] 2) Soxhlet Extraction Method ^[2]
16	pH	Electrometric method ^[2]
17	Phenols	Distillation, Direct Photometric Method ^[2]
18	Selenium	Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[2]
19	Sulfide	1) Iodometric Method ^[2] 2) Methylene Blue Method ^[2]
20	Temperature	Laboratory and Field Methods ^[2]
21	Total Dissolved Solids	Dried at 180 °C ^[2]
22	Total Kjeldahl Nitrogen	Macro-Kjeldahl Method ^[2]
23	Total Suspended Solids	Dried at 103-105 °C ^[2]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
24	Trivalent Chromium	Digestion, Direct Air-Acetylene Flame Method; Colorimetric Method; Calculation ^[2]
25	Zinc	Digestion, Direct Air-Acetylene Flame Method ^[2]

เอกสารอ้างอิง

1. สมาคมวิศวกรรมสิ่งแวดล้อมแห่งประเทศไทย. คู่มือวิเคราะห์น้ำเสีย. พิมพ์ครั้งที่ 4. กรุงเทพฯ: เรือนแก้วการพิมพ์, 2547
2. APHA, AWWA, WEF. Standard Methods for the Examination of Water and Wastewater. 23rd ed. Washington, DC: APHA, 2017.

คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน

วันที่ 10 เดือน เมษายน พ.ศ. 2567

เรียน อธิบดีกรมโรงงานอุตสาหกรรม

ข้าพเจ้า ☐ ผู้รับใบอนุญาตประกอบกิจการโรงงาน.....

☒ บริษัท/ห้างหุ้นส่วนจำกัด /.....บริษัท เอ็นไวรอนเม้นท์ แอนด์ แกลบอราตอรี.....

สถาน.....

ได้รับทราบ อธิบดีกรมโรงงานอุตสาหกรรม เรื่อง ห้องปฏิบัติการวิเคราะห์เอกชน พ.ศ. 2566 เคยตลอดแล้ว และยินยอมปฏิบัติตาม
ประกาศฯ ทุกประการ และได้แนบเอกสารต่างๆ ตามรายการเอกสารประกอบการพิจารณา (แบบ ปอ.1-1) มาพร้อมนี้

รายการขอดำเนินการ

การดำเนินการ	จำนวนสารมลพิษ					
	น้ำเสีย/น้ำทิ้ง (รายการ)	น้ำใต้ดิน (รายการ)	อากาศ (รายการ)	สิ่งปฏิกูลหรือ วัสดุที่ไม่ใช้ แล้ว(รายการ)	ดิน (รายการ)	รวมทั้งสิ้น (รายการ)
<input type="checkbox"/> ขอขึ้นทะเบียนห้องปฏิบัติการ วิเคราะห์เอกชน						
<input checked="" type="checkbox"/> ต่ออายุห้องปฏิบัติการ วิเคราะห์เอกชน	26					
<input type="checkbox"/> เปลี่ยนแปลงสารมลพิษที่ วิเคราะห์ ○ เพิ่มสารมลพิษ ○ ยกเลิกสารมลพิษ						
<input type="checkbox"/> เปลี่ยนแปลงบุคลากร ○ เพิ่มบุคลากร ○ ยกเลิกบุคลากร	จำนวน.....ราย (รายละเอียดตาม แบบ ปว.1) จำนวน.....ราย (รายละเอียดตาม แบบ ปว.1-1)					
<input type="checkbox"/> ยกเลิกห้องปฏิบัติการวิเคราะห์เอกชน						
<input type="checkbox"/> อื่นๆ โปรดระบุ.....						



ที่ อก ๐๓๒๐/ [REDACTED]

กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๕๐๐

๑ พฤศจิกายน ๒๕๖๕

เรื่อง ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

เรียน กรรมการผู้จัดการ บริษัท เอสจีเอส (ประเทศไทย) จำกัด (สาขาระยอง)

อ้างถึง คำขอต่ออายุของห้องปฏิบัติการวิเคราะห์เอกชน ลงวันที่ ๑๒ กันยายน ๒๕๖๕

สิ่งที่ส่งมาด้วย เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
บริษัท เอสจีเอส (ประเทศไทย) จำกัด (สาขาระยอง) จำนวน ๒๒ แผ่น

ตามหนังสือที่อ้างถึง บริษัท เอสจีเอส (ประเทศไทย) จำกัด (สาขาระยอง) ขอต่ออายุหนังสือ
รับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ว-๑๙๗ สถานที่ตั้งเลขที่ ๑/๒๐๙ และ ๑/๒๑๑ หมู่ที่ ๑
ตำบลบ้านฉาง อำเภอบ้านฉาง จังหวัดระยอง ต่อกรมโรงงานอุตสาหกรรม นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้บริษัท เอสจีเอส (ประเทศไทย) จำกัด (สาขาระยอง)
ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน โดยมีองค์ประกอบดังนี้






ค. ขอบขายสารมลพิษที่เตรียมทะเบียนเหวเคราะห์เนนาเสีย จำนวน ๔๔ รายการ
น้ำใต้ดิน จำนวน ๑๒๓ รายการ อากาศเสีย (ปล่อยระบาย) จำนวน ๒๘ รายการ สิ่งปฏิภูลหรือวัสดุที่ไม่ใช่แล้ว
จำนวน ๓๗ รายการ และดิน จำนวน ๑๒๓ รายการ รวมทั้งสิ้นจำนวน ๓๕๕ รายการ ตามสิ่งที่ส่งมาด้วย

หนังสือฉบับนี้จะหมดอายุในวันที่ ๑๒ ตุลาคม ๒๕๖๘ หากประสงค์จะต่ออายุหนังสือ
รับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ให้ยื่นคำขอต่ออายุพร้อมเอกสารประกอบคำขอต่อกรมโรงงาน
อุตสาหกรรมภายใน ๓๐ วัน ก่อนวันสิ้นอายุของหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ซึ่งคำขอ
ต่ออายุดังกล่าวขอรับได้ที่กรมโรงงานอุตสาหกรรม ทั้งนี้ สามารถยื่นคำขอผ่านระบบอิเล็กทรอนิกส์ได้ที่หน้า
เว็บไซต์กรมโรงงานอุตสาหกรรม ตาม QR Code ท้ายหนังสือนี้

จึงเรียนมาเพื่อทราบ



ศูนย์วิจัยและเตือนภัยมลพิษโรงงานภาคตะวันออก

โทร. ๐ ๓๓๑๓ ๖๐๕๙ ต่อ ๕๐๐๑-๒

ไปรษณีย์อิเล็กทรอนิกส์ eirw@diw.mail.go.th



ยื่นคำขอผ่านระบบอิเล็กทรอนิกส์



“อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว”



เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท เอสจีเอส (ประเทศไทย) จำกัด (สาธารณะยง)

เลขทะเบียน ว-๑๙๗

ที่ อก ๐๓๒๐/๑๖๐๔๑

ลงวันที่ ๑ พฤศจิกายน ๒๕๖๕

ขอข่ายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๓๕๕ รายการ
น้ำเสีย จำนวน 44 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Aldrin	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
2	Arsenic	Digestion, Inductively Coupled Plasma Method ^[4]
3	Barium	Digestion, Inductively Coupled Plasma Method ^[4]
4	α -BHC	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
5	β -BHC	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
6	δ -BHC	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
7	γ -BHC	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
8	Biochemical Oxygen Demand	5-Day BOD Test, Membrane Electrode Method ^[4]
9	Cadmium	Digestion, Inductively Coupled Plasma Method ^[4]
10	Chemical Oxygen Demand	Closed Reflux, Titrimetric Method ^[4]
11	Chlordane	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
12	Color	ADMI Weighted – Ordinate Spectrophotometric Method ^[4]
13	Copper	Digestion, Inductively Coupled Plasma Method ^[4]
14	Cyanide	Distillation, Colorimetric Method ^[4]
15	p,p'-DDD	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
16	p,p'-DDE	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
17	o,p'-DDT	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
18	p,p'-DDT	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
19	Dieldrin	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
20	Endosulfan I	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
21	Endosulfan II	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
22	Endosulfan Sulfate	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
23	Endrin	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
24	Endrin Aldehyde	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
25	Formaldehyde	Distillation, Colorimetric Method ^[3]
26	Heptachlor	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
27	Heptachlor Epoxide	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
28	Hexavalent Chromium	Filtration, Colorimetric Method ^[4]
29	Lead	Digestion, Inductively Coupled Plasma Method ^[4]
30	Manganese	Digestion, Inductively Coupled Plasma Method ^[4]
31	Mercury	Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[4]
32	Methoxychlor	Liquid-Liquid Extraction, Gas Chromatographic Method ^[4]
33	Nickle	Digestion, Inductively Coupled Plasma Method ^[4]
34	Oil and Grease	Liquid-Liquid, Partition-Gravimetric Method ^[4]
35	pH	Electrometric Method ^[4]
36	Phenols	Distillation, Direct Photometric Method ^[4]
37	Selenium	Digestion, Inductively Coupled Plasma Method ^[4]
38	Temperature	Field Method ^[4]
39	Total Chromium	Digestion, Inductively Coupled Plasma Method ^[4]
40	Total Dissolved Solids	Dried at 180 °C ^[4]
41	Total Kjeldahl Nitrogen	Digestion, Distillation, Titrimetric Method ^[4]
42	Total Suspended Solids	Dried at 103-105 °C ^[4]
43	Trivalent Chromium	Digestion, Inductively Coupled Plasma Method; Filtration, Colorimetric Method. Calculation ^[4]
44	Zinc	Digestion, Inductively

น้ำใต้ดิน จำนวน 123 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Acenaphthene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
2	Acetone	Purge and Trap Gas Chromatographic / Mass Spectrometric Method ^[4]
3	Aldrin	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
4	Anthracene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
5	Antimony	Digestion, Inductively Coupled Plasma Method ^[4]
6	Arsenic	Digestion, Inductively Coupled Plasma Method ^[4]
7	Atrazine	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
8	Barium	Digestion, Inductively Coupled Plasma Method ^[4]
9	Benzene	Purge and Trap Gas Chromatographic / Mass Spectrometric Method ^[4]
10	Benzo(a)anthracene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
11	Benzo(b)fluoranthene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
12	Benzo(k)fluoranthene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
13	Benzoic acid	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
14	Benzo(a)pyrene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
15	Benzo(g,h,i)perylene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
16	Beryllium	Digestion, Inductively Coupled Plasma Method ^[4]
17	Bis(2-chloroethyl)ether	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
18	Bis(2-Ethylhexyl)phthalate	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
19	Bromodichloromethane	Purge and Trap Gas Chromatographic / Mass Spectrometric Method ^[4]
20	Bromoform	Purge and Trap, Gas Chromatographic / Mass Spectrometric Meth

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
21	Butyl benzyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
22	Cadmium	Digestion, Inductively Coupled Plasma Method ^[4]
23	Carbazole	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
24	Carbon disulfide	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
25	Carbon tetrachloride	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
26	Chlordane	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
27	p-Chloroaniline	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
28	Chlorobenzene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
29	Chlorodibromomethane	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
30	Chloroform	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
31	2-Chlorophenol	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
32	Chromium	Digestion, Inductively Coupled Plasma Method ^[4]
33	Chromium Hexavalent	Filtration, Colorimetric Method ^[4]
34	Chromium Trivalent	Digestion, Inductively Coupled Plasma Method ; Filtration, Colorimetric Method; Calculation ^[4]
35	Chrysene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
36	Cyanide	Distillation, Colorimetric Method
37	2,4-D	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
38	DDD	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
39	DDE	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
40	DTT	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
41	Dibenz(a,h)anthracene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
42	Di-n-Butyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
43	1,2-Dichlorobenzene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
44	1,3-Dichlorobenzene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
45	1,4-Dichlorobenzene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
46	3,3-Dichlorobenzidine	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
47	1,1-Dichloroethane	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
48	1,2-Dichloroethane	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
49	1,1-Dichloroethylene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
50	cis-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
51	trans-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
52	2,4-Dichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
53	1,2-Dichloropropane	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
54	1,3-Dichloropropane	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
55	1,3-Dichloropropene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
56	Dieldrin	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
57	Diethyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
58	2,4-Dimethylphenol	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]

2,4-Dinitrophenol...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
59	2,4-Dinitrophenol	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
60	2,4-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
61	2,6-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
62	Di-n-octyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
63	Endosulfan	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
64	Endrin	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
65	Ethylbenzene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
66	Fluoranthene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
67	Fluorene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
68	Heptachlor	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
69	Heptachlor epoxide	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
70	Hexachlorobenzene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
71	Hexachloro-1,3-butadiene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
72	α -HCH	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
73	β -HCH	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
74	γ -HCH	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
75	Hexachlorocyclopentadiene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
76	Hexachloroethane	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
77	n-Hexane	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method
78	Indeno(1,2,3-cd)pyrene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
79	Isophorone	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
80	Lead	Digestion, Inductively Coupled Plasma Method ^[4]
81	Manganese	Digestion, Inductively Coupled Plasma Method ^[4]
82	Mercury	Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[4]
83	Methoxychlor	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
84	Methyl Bromide	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
85	Methylene Chloride	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
86	2-Methylnaphthalene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
87	2-Methylphenol	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
88	Methyl tert-butyl ether	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
89	Naphthalene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
90	Nickel	Digestion, Inductively Coupled Plasma Method ^[4]
91	Nitrobenzene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
92	N-Nitrosodiphenylamine	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
93	N-Nitrosodi-n-propylamine	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
94	Pentachlorophenol	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
95	pH	Electrometric Method ^[4]
96	Phenanthrene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
97	Phenol	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
98	Pyrene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
99	Selenium	Digestion, Inductively Coupled Plasma Method ^[4]
100	Silver	Digestion, Inductively Coupled Plasma Method ^[4]
101	Styrene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
102	1,1,2,2-Tetrachloroethane	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
103	Tetrachloroethylene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
104	Toluene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
105	Toxaphene	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
106	TPH (C ₅ -C ₈)	Purge and Trap, Gas Chromatographic Mass Spectrometric Method
107	TPH (C ₈ -C ₁₆)	Purge and Trap, Gas Chromatographic Mass Spectrometric Method
108	TPH (C ₁₆ -C ₃₅)	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
109	1,2,4-Trichlorobenzene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
110	1,1,1-Trichloroethane	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
111	1,1,2-Trichloroethane	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
112	Trichloroethylene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
113	2,4,5-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
114	2,4,6-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic / Mass Spectrometric Method ^[4]
115	1,3,5-Trimethylbenzene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
116	Vanadium	Digestion, Inductively Coupled Plasma Method ^[4]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
117	Vinyl acetate	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
118	Vinyl chloride	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
119	m-Xylene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
120	o-Xylene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
121	p-Xylene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
122	Xylene (Total)	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[4]
123	Zinc	Digestion, Inductively Coupled Plasma Method ^[4]

อากาศเสีย (ปล่องระบาย) จำนวน 28 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Antimony	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
2	Arsenic	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
3	Beryllium	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
4	Cadmium	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
5	Carbon Monoxide	Instrumental Analyzer Method ^[5]
6	Chlorine	Isokinetic Sampling, Ion Chromatographic Method ^[7]
7	Chromium	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
8	Cobalt	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
9	Copper	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
10	Cresol	Absorption Sampling, Gas Chromatographic Method ^[6]
11	Dioxin/Furans	Isokinetic Sampling, Analysis by ISO/IEC 17025 Accredited Laboratory

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
12	Hydrogen Chloride	Isokinetic Sampling, Ion Chromatographic Method ^[7]
13	Hydrogen Fluoride	Isokinetic Sampling, Ion Chromatographic Method ^[7]
14	Hydrogen Sulfide	Absorption Sampling, Iodometric Method ^[5]
15	Lead	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
16	Manganese	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
17	Mercury	Isokinetic Sampling, Digestion, Cold-Vapour Atomic Absorption Spectrometric Method ^[5]
18	Nickel	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
19	Opacity	Ringelmann's Method ^[1]
20	Oxides of Nitrogen	1) Absorption Sampling, Colorimetric Method ^[6] 2) Instrumental Analyzer Method ^[7]
21	Tellurium	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
22	Tin	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
23	Total Suspended Particulate	Isokinetic Sampling, Gravimetric Method ^[7]
24	Selenium	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
25	Sulfur Dioxide	1) Absorption Sampling, Barium-Thorin Titrimetric Method ^[5] 2) Instrumental Analyzer Method ^[5]
26	Sulfuric Acid	Isokinetic Sampling, Barium-Thorin Titrimetric Method ^[6]
27	Vanadium	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
28	Xylene	Adsorption Sampling, Gas Chromatographic Method ^[6]

สิ่งปฏิกูลหรือวัสดุที่ไม่ใช้แล้ว จำนวน 37 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Aldrin	Ultrasonic Extraction, Gas Chromatographic Method ^[12,20,21]
2	Antimony	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[2,15] 2) Digestion, Inductively Coupled Plasma Method ^[8,15]
3	Arsenic	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[2,15] 2) Digestion, Inductively Coupled Plasma Method ^[8,15]
4	Barium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[2,15] 2) Digestion, Inductively Coupled Plasma Method ^[8,15]
5	Beryllium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[2,15] 2) Digestion, Inductively Coupled Plasma Method ^[8,15]
6	Cadmium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[2,15] 2) Digestion, Inductively Coupled Plasma Method ^[8,15]
7	Chlordane	Ultrasonic Extraction, Gas Chromatographic Method ^[12,20,21]
8	Chromium (III)	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method; Waste Extraction Colorimetric Method; Calculation ^[10,17] 2) Digestion, Inductively Coupled Plasma Method; Alkaline Digestion, Colorimetric Method; Calculation Method ^[2,10,17]
9	Chromium (VI)	1) Waste Extraction, Digestion, Colorimetric Method ^[10,17] 2) Alkaline Digestion, Colorimetric Method ^[10,17]
10	Cobalt	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[2,15] 2) Digestion, Inductively Coupled Plasma Method ^[8,15]
11	Copper	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[2,15] 2) Digestion, Inductively Coupled Plasma Method ^[8,15]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
12	Dieldrin	Ultrasonic Extraction, Gas Chromatographic Method ^[12,20,21]
13	DDD	Ultrasonic Extraction, Gas Chromatographic Method ^[12,20,21]
14	DDE	Ultrasonic Extraction, Gas Chromatographic Method ^[12,20,21]
15	DDT	Ultrasonic Extraction, Gas Chromatographic Method ^[12,20,21]
16	2,4-D (2,4-Dichlorophenoxyacetic acid)	Ultrasonic Extraction, Gas Chromatographic Method ^[12,20,21]
17	Endrin	Ultrasonic Extraction, Gas Chromatographic Method ^[12,20,21]
18	Heptachlor	Ultrasonic Extraction, Gas Chromatographic Method ^[12,20,21]
19	Kepone	Ultrasonic Extraction, Gas Chromatographic Method ^[12,20,21]
20	Lead	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[2,11] 2) Digestion, Inductively Coupled Plasma Method ^[8,15]
21	Lindane	Ultrasonic Extraction, Gas Chromatographic Method ^[12,20,21]
22	Mercury	1) Waste Extraction, Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[2,18] 2) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[18]
23	Methoxychlor	Ultrasonic Extraction, Gas Chromatographic Method ^[12,20,21]
24	Mirex	Ultrasonic Extraction, Gas Chromatographic Method ^[12,20,21]
25	Molybdenum	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[2,15] 2) Digestion, Inductively Coupled Plasma Method ^[8,15]
26	Polychlorinated Biphenyls (PCBs)	Ultrasonic Extraction, Gas Chromatographic Method ^[12,20,21]
27	Pentachlorophenol	Ultrasonic Extraction, Gas Chromatographic Method ^[12,20,21]
28	Nickel	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[2,15] 2) Digestion, Inductively Coupled Plasma Method ^[8,15]
29	Selenium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[2,15] 2) Digestion, Inductively Coupled Plasma Method ^[8,15]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
30	Silver	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[2,15] 2) Digestion, Inductively Coupled Plasma Method ^[8,15]
31	Silvex; 2,4,5-Trichlorophenoxypropionic acid	Ultrasonic Extraction, Gas Chromatographic Method ^[12,20,21]
32	Thallium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[2,15] 2) Digestion, Inductively Coupled Plasma Method ^[8,15]
33	Total Chromium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method; Waste Extraction Colorimetric Method; Calculation ^[10,17] 2) Digestion, Inductively Coupled Plasma – Atomic Emission Spectrometry Method ^[8,15]
34	Toxaphene	Ultrasonic Extraction, Gas Chromatographic Method ^[12,20,21]
35	Trichloroethylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[14,22]
36	Vanadium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[2,15] 2) Digestion, Inductively Coupled Plasma Method ^[8,15]
37	Zinc	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[2,15] 2) Digestion, Inductively Coupled Plasma Method ^[8,15]

ดิน จำนวน 123 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Acenaphthene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
2	Acetone	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
3	Aldrin	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
4	Anthracene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
5	Antimony	Digestion, Inductively Coupled Plasma Method ^[9,15]

6 Arsenic...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
6	Arsenic	Digestion, Inductively Coupled Plasma Method ^[9,15]
7	Atrazine	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
8	Barium	Digestion, Inductively Coupled Plasma Method ^[9,15]
9	Benzo(a)anthracene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
10	Benzene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
11	Benzo(b)fluoranthene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
12	Benzo(k)fluoranthene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
13	Benzoic acid	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
14	Benzo(a)pyrene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
15	Benzo(g,h,i)perylene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
16	Beryllium	Digestion, Inductively Coupled Plasma Method ^[9,15]
17	Bis(2-Chloroethyl)ether	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
18	Bis(2-Ethylhexyl)phthalate	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
19	Bromodichloromethane	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
20	Bromoform	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
21	Butyl benzyl phthalate	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
22	Cadmium	Digestion, Inductively Coupled Plasma Method ^[9,15]
23	Carbazole	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
24	Carbon disulfide	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
25	Carbon tetrachloride	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
26	Chlordane	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
27	p-Chloroaniline	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
28	Chlorobenzene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
29	Chlorodibromomethane	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
30	Chloroform	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
31	2-Chlorophenol	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
32	Chromium	Digestion, Inductively Coupled Plasma Method ^[9,15]
33	Chromium (III)	Digestion, Inductively Coupled Plasma Method ; Filtration, Colorimetric Method; Calculation ^[9,10,15]
34	Chromium (VI)	Alkaline Digestion, Colorimetric Method ^[10]
35	Chrysene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
36	Cyanide	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
37	2,4-D	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
38	DDD	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
39	DDE	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
40	DDT	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
41	Dibenz(a,h)anthracene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
42	Di-n-Butyl phthalate	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
43	1,2-Dichlorobenzene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
44	1,3-Dichlorobenzene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
45	1,4-Dichlorobenzene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
46	3,3-Dichlorobenzidine	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
47	1,1-Dichloroethane	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
48	1,2-Dichloroethane	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
49	1,1-Dichloroethylene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
50	cis-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
51	trans-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
52	2,4-Dichlorophenol	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
53	1,2-Dichloropropane	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
54	1,3-Dichloropropane	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
55	1,3-Dichloropropene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
56	Dieldrin	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
57	Diethyl phthalate	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
58	2,4-Dimethylphenol	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
59	2,4-Dinitrophenol	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
60	2,4-Dinitrotoluene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
61	2,6-Dinitrotoluene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
62	Di-n-octyl phthalate	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
63	Endosulfan	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
64	Endrin	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
65	Ethylbenzene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
66	Fluoranthene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
67	Fluorene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
68	Heptachlor	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
69	Heptachlor epoxide	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
70	Hexachlorobenzene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
71	Hexachloro-1,3-butadiene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
72	α -HCH	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
73	β -HCH	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
74	γ -HCH	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
75	Hexachlorocyclopentadiene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
76	Hexachloroethane	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
77	n-Hexane	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
78	Indeno(1,2,3-cd)pyrene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
79	Isophorone	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
80	Lead	Digestion, Inductively Coupled Plasma Method ^[9,15]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
81	Manganese	Digestion, Inductively Coupled Plasma Method ^[9,15]
82	Mercury	Digestion, Cold vapor Atomic Absorption Spectrometric Method
83	Methoxychlor	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
84	Methyl Bromide	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
85	Methylene Chloride	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
86	2-Methylnaphthalene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
87	2-Methylphenol	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
88	Methyl tert-butyl ether	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
89	Naphthalene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
90	Nickel	Digestion, Inductively Coupled Plasma Method ^[9,15]
91	Nitrobenzene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
92	N-Nitrosodiphenylamine	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
93	N-Nitrosodi-n-propylamine	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
94	Pentachlorophenol	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
95	Phenanthrene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
96	Phenol	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
97	Polychlorinated Biphenyls (PCBs)	Ultrasonic Extraction, Gas Chromatographic Method ^[9,16,17]
98	Pyrene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[23,24]
99	Selenium	Digestion, Inductively Coupled Plasma Method ^[9,15]
100	Silver	Digestion, Inductively Coupled Plasma Method ^[9,15]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
101	Styrene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
102	1,1,2,2-Tetrachloroethane	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
103	Tetrachloroethylene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
104	Toluene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
105	Toxaphene	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[9,10]
106	TPH (C ₅ -C ₈)	Purge and Trap, Gas Chromatographic Mass Spectrometric Method ^[14,22]
107	TPH (C ₈ -C ₁₆)	Ultrasonic Extraction, Gas Chromatographic Mass Spectrometric Method ^[9,10,18]
108	TPH (C ₁₆ -C ₃₅)	Ultrasonic Extraction, Gas Chromatographic Mass Spectrometric Method ^[10,18]
109	1,2,4-Trichlorobenzene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
110	1,1,1-Trichloroethane	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
111	1,1,2-Trichloroethane	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
112	Trichloroethylene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
113	2,4,5-Trichlorophenol	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[9,10]
114	2,4,6-Trichlorophenol	Microwave Extraction, Gas Chromatographic / Mass Spectrometric Method ^[9,10]
115	1,3,5-Trimethylbenzene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
116	Vanadium	Digestion, Inductively Coupled Plasma Method ^[9,15]
117	Vinyl Acetate	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[14,22]
118	Vinyl Chloride	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
119	m-Xylene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[6,8]
120	o-Xylene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[6,8]
121	p-Xylene	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[6,8]
122	Xylene (Total)	Purge and Trap, Gas Chromatographic / Mass Spectrometric Method ^[6,8]
123	Zinc	Digestion, Inductively Coupled Plasma Method ^[9,15]

เอกสารอ้างอิง


1. กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม พ.ศ.2549 เรื่อง กำหนดค่าปริมาณ
เขม่าควันที่เจือปนในอากาศที่ระบายออกจากปล่องของหม้อน้ำโรงสีข้าวที่ใช้กลบเป็นเชื้อเพลิง.
ราชกิจจานุเบกษา. 4 ธันวาคม 2549. เล่มที่ 123 ตอนพิเศษ 125 ง.
2. กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม พ.ศ.2548 เรื่อง การกำจัดสิ่งปฏิกูล
หรือวัสดุที่ไม่ใช้แล้ว. ราชกิจจานุเบกษา. 25 มกราคม 2549. เล่มที่ 123 ตอนพิเศษ 11ง.
3. สมาคมวิศวกรสิ่งแวดล้อมแห่งประเทศไทย. คู่มือวิเคราะห์น้ำเสีย. พิมพ์ครั้งที่ 4. กรุงเทพฯ:
เรือนแก้วการพิมพ์, 2547
4. APHA, AWWA, WEF. **Standard Methods for the Examination of Water and
Wastewater.** 23rd ed. Washington, DC : APHA, 2017
5. United States Environmental Protection Agency. **Standards of Performance
for New Stationary Sources.** 40 CFR 60. Appendix A, 2017
6. United States Environmental Protection Agency. **Standards of Performance
for New Stationary Sources.** 40 CFR 60. Appendix A, 2019
7. United States Environmental Protection Agency. **Standards of Performance
for New Stationary Sources.** 40 CFR 60. Appendix A, 2020
8. United States Environmental Protection Agency. **Test Methods for Evaluation
Solid Waste Physical/Chemical Methods. Acid Digestion of Sediments Sludges, and Soils.
SW-846 Method 3050B,** 1996.
9. United States Environmental Protection Agency. **Test Methods for Evaluation
Solid Waste Physical/Chemical Methods. Acid Digestion of Sediments, Sludges, and Soils.
SW-846 Method 3051A,** 2007
10. United States Environmental Protection Agency. **Test Methods for
Evaluation Solid Waste Physical/Chemical Methods. Alkaline Digestion for Hexavalent
Chromium. SW-846 Method 3060A,** 1996.

11. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. SW-846, 2006.
12. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Ultrasonic Extraction. SW-846 Method 3550C, 2007.
13. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Closed-System Purge-and-Trap and Extraction for Volatile Organics in Soil and Waste Samples. SW-846 Method 5035A, 2002
14. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Closed-System Purge-and-Trap for Aqueous Samples. SW-846 Method 5035C, 2003.
15. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Inductively Coupled Plasma – optical Emission Spectrometry. SW-846 Method 6010D, 2018
16. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Antimony and Arsenic (Atomic Absorption, Borohydride Reduction). SW-846 Method 7062A, 1994.
17. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Chromium, Hexavalent (Colorimetric). Method 7196A, 1992.
18. United States Environment Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Mercury in Solid or Semisolid Waste (Manual Cold-Vapor Technique). SW-846 Method 7471B, 2007.
19. United States Environment Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Nonhalogenated Organics Using GC/FID. SW-846 Method 8015D, 2003.
20. United States Environment Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Organochlorine Pesticide by Gas Chromatography. SW-846 Method 8081B, 2007.
21. United States Environment Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Polychlorinated Biphenyls (PCBs) By Gas Chromatography. SW-846 Method 8082A, 2007
22. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Volatile Organic Compounds by Gas Chromatography/Mass Spectrometry (GC/MS). SW-846 Method 8260D, 2018.

23. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Semivolatile Organic Compounds by Gas Chromatography/Mass Spectrometry. SW-846 Method 8270E, 2018.

24. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Microwave Extraction, Gas Chromatography/Mass Spectrometry. SW-846 Method 3546, 2007.

25. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Cyanide Extraction Procedure for Solids and Oils. SW-846 Method 9013A, 2014.



ที่ อก ๐๓๑๐(๑)/



กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๖๙ มกราคม ๒๕๖๗

เรื่อง ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

เรียน กรรมการผู้จัดการ บริษัท บุโร เวิร์ทส์ เอคิว แล็บ (ประเทศไทย) จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๑๖ ตุลาคม ๒๕๖๖

สิ่งที่ส่งมาด้วย เอกสารแนบท้ายหนังสือต่ออายุรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
บริษัท บุโร เวิร์ทส์ เอคิว แล็บ (ประเทศไทย) จำกัด จำนวน ๒ แผ่น

ตามหนังสือที่อ้างถึง บริษัท บุโร เวิร์ทส์ เอคิว แล็บ (ประเทศไทย) จำกัด ห้องปฏิบัติการ
วิเคราะห์เอกชน เลขทะเบียน ว-๓๑๓ สถานที่ตั้งเลขที่ ๑๑๑ หมู่ที่ ๙ อาคารสำนักงานกลาง อุทยานวิทยาศาสตร์
แห่งประเทศไทย ชั้นที่ ๑ และชั้นที่ ๒ ห้องเลขที่ P-๑๐๓ และ P-๒๐๔C-E ถนนพหลโยธิน ตำบลคลองหนึ่ง
อำเภอคลองหลวง จังหวัดปทุมธานี ต่อกรมโรงงานอุตสาหกรรม นั้น

ต่ออายุหนังสือ

รับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ให้ยื่นคำขอต่ออายุพร้อมเอกสารประกอบคำขอต่อกรมโรงงาน
อุตสาหกรรมภายใน ๓๐ วัน ก่อนวันสิ้นอายุของหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ทั้งนี้
สามารถยื่นคำขอผ่านระบบอิเล็กทรอนิกส์ได้ที่หน้าเว็บไซต์กรมโรงงานอุตสาหกรรม

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

กองวิจัยและเตือนภัยมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕

โทรสาร ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๙๙

ไปรษณีย์อิเล็กทรอนิกส์ saraban@diw.mail.go.th



“อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว”



เอกสารแนบท้ายหนังสือต่ออายุรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท บุโร เวอร์ทีส เอคิว แล็บ (ประเทศไทย) จำกัด เลขทะเบียน ว-๓๑๓

ที่ ออก ๐๓๑๐(๑)/ [REDACTED] ลงวันที่ ๑๙ มกราคม ๒๕๖๗

ขอขยายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๓๕ รายการ

น้ำเสีย จำนวน 35 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Aldrin	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method
2	Arsenic	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method
3	Barium	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method
4	α -BHC	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method
5	β -BHC	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method
6	δ -BHC	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method
7	γ -BHC	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method
8	Cadmium	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method
9	Chlordane	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method
10	Chromium	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method
11	Color	ADMI Weighted-Ordinate Spectrophotometric Method
12	Copper	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method
13	o,p'-DDT	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method
14	4,4'-DDD	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method
15	4,4'-DDE	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method
16	4,4'-DDT	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method
17	Dieldrin	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method [REDACTED]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
18	Endosulfan I	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method
19	Endosulfan II	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method
20	Endosulfan sulfate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method
21	Endrin	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method
22	Heptachlor	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method
23	Heptachlor epoxide	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method
24	Lead	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method
25	Manganese	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method
26	Mercury	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method
27	Mirex	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method
28	Methoxychlor	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method
29	Nickel	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method
30	pH	Electrometric Method
31	Selenium	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method
32	Total Dissolved Solids	Dried at 180 °C
33	Total Kjeldahl Nitrogen	Macro-Kjeldahl Method
34	Total Suspended Solids	Dried from 103 to 105 °C
35	Zinc	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method

เอกสารอ้างอิง

APHA, AWWA, WEF. Standard Methods for the Examination of Water and Wastewater. 24th ed. Washington, DC: APHA, 2023.