

ภาคผนวก ณ

ใบรับรองการสอบเทียบเครื่องมือ



right solutions.  
right partner.

รายการเครื่องมือที่ใช้ในการวิเคราะห์ / ทดสอบ

Sample Name	Parameter	Equipment Name	ID No.	Calibrated Date	Next Cal	Freq. Calibrate (Months)
Noise	Leq 24 hrs	Sound Calibrator	BKK_FS0631	26-Jan-24	25-Jan-25	18
Noise	Leq 24 hrs	Sound Level Meter	BKK_FS0034	29-Nov-23	28-Nov-24	12
Noise	Leq 24 hrs	Sound Level Meter	BKK_FS0032	1-Sep-23	1-Sep-24	12
Noise	Leq 24 hrs	Sound Level Meter	BKK_FS0098	9-Jul-24	9-Jul-25	12
Noise	Leq 24 hrs	Sound Calibrator	BKK_FS0030	27-Jan-24	27-Jan-25	12
Noise	Leq 24 hrs	Sound Level Meter	BKK_FS0099	30-Aug-24	30-Aug-25	12
Noise	Leq 24 hrs	Sound Level Meter	BKK_FS0100	30-Aug-24	30-Aug-25	12
Noise	Leq 24 hrs	Sound Level Meter	BKK_FS0101	29-May-24	29-May-25	12
Water Lab	pH at 25 °C	pH meter	BKK_EN0342	17-Oct-24	17-Oct-25	12
Water Lab	Total Hardness	Burette	BKK_EN0171	27-Feb-24	27-Aug-25	18
Water Lab	Total Hardness	Burette	BKK_EN0171	27-Feb-24	27-Aug-25	18
Water Lab	Color	Chamber (Cold Room)	BKK_EN0167	6-Dec-23	6-Jun-25	18
Water Lab	Turbidity	Chamber (Cold Room)	BKK_EN0167	6-Dec-23	6-Jun-25	18
Water Lab	Dissolved Oxygen	Burette	BKK_EN0171	27-Feb-24	27-Aug-25	18
Water Lab	Dissolved Oxygen	Chamber (Cold Room)	BKK_EN0167	6-Dec-23	6-Jun-25	18
Water Lab	Nitrate	Ion Chromatography	BKK_EN0069	12-Jan-24	12-Jan-25	12
Water Lab	Phosphate	Ion Chromatography	BKK_EN0069	12-Jan-24	12-Jan-25	12
Water Lab	Chloride	Ion Chromatography	BKK_EN0069	12-Jan-24	12-Jan-25	12
Water Lab	Sulfate	Ion Chromatography	BKK_EN0069	12-Jan-24	12-Jan-25	12
Water Lab	Total Suspended Solids	Electronic Top-Loading Balance	BKK_EN0003	2-Aug-24	2-Aug-25	12
Water Lab	Total Suspended Solids	Oven	BKK_EN0273	14-May-24	14-Nov-25	18
Water Lab	Total Solids	Electronic Top-Loading Balance	BKK_EN0003	2-Aug-24	2-Aug-25	12
Water Lab	Total Solids	Oven	BKK_EN0273	14-May-24	14-Nov-25	18
Water Lab	Total Dissolved Solids 180°C	Electronic Top-Loading Balance	BKK_EN0003	2-Aug-24	2-Aug-25	12
Water Lab	Total Dissolved Solids 180°C	Oven	BKK_EN0273	14-May-24	14-Nov-25	18
Water Lab	Total Kjeldahl Nitrogen	Digestion Unit	BKK_EN0366	21-Apr-24	21-Apr-25	12
Water Lab	Total Kjeldahl Nitrogen	Discrete analyzer	BKK_EN0037	16-Aug-24	16-Aug-25	12
Water Lab	Conductivity	Conductivity meter	BKK_EN0373	25-Dec-23	25-Dec-24	12
Water Lab	Salinity	Conductivity meter	BKK_EN0373	25-Dec-23	25-Dec-24	12
Water Lab	BOD	DO Meter	BKK_EN0017	16-Nov-23	16-May-25	18
Water Lab	BOD	Incubator	BKK_EN0304	20-Mar-24	20-Mar-25	12
Water Lab	BOD	Burette	BKK_EN0171	27-Feb-24	27-Aug-25	18
Water Lab	COD	Hot Block	BKK_EN0222	22-Apr-24	22-Apr-25	12
Water Lab	COD	Spectrophotometer	BKK_EN0018	13-Sep-24	13-Sep-25	12
Water Lab	Oil & Grease	Electronic Top-Loading Balance	BKK_EN0003	2-Aug-24	2-Aug-25	12
Water Lab	Oil & Grease	Water Bath	BKK_EN0148	6-Jul-23	4-Jan-25	18
Water Lab	Temperature	Digital Thermometer With Sensor	BKK_LG0056	22-Dec-23	22-Dec-24	12
Water Lab	Hexavalent Chromium	Spectrophotometer	BKK_EN0018	13-Sep-24	13-Sep-25	12
Water Lab	Iron	ICP-AES	BKK_EL0043	8-Oct-24	3-Apr-26	18
Water Lab	Iron	Hot Block	BKK_EL0054	22-Sep-23	22-Mar-25	18
Water Lab	Iron	Chamber (Cold Room)	BKK_EN0167	6-Dec-23	6-Jun-25	18
Water Lab	Lead	ICP-AES	BKK_EL0043	8-Oct-24	3-Apr-26	18
Water Lab	Lead	Hot Block	BKK_EL0054	22-Sep-23	22-Mar-25	18
Water Lab	Lead	Chamber (Cold Room)	BKK_EN0167	6-Dec-23	6-Jun-25	18
Water Lab	Manganese	ICP-AES	BKK_EL0043	8-Oct-24	3-Apr-26	18
Water Lab	Manganese	Hot Block	BKK_EL0054	22-Sep-23	22-Mar-25	18
Water Lab	Manganese	Chamber (Cold Room)	BKK_EN0167	6-Dec-23	6-Jun-25	18
Water Lab	Calcium	ICP-OES	BKK_EL0037	29-Feb-24	28-Feb-25	12
Water Lab	Calcium	Hot Block	BKK_EL0054	22-Sep-23	22-Mar-25	18
Water Lab	Calcium	Chamber (Cold Room)	BKK_EN0167	6-Dec-23	6-Jun-25	18
Water Lab	Magnesium	ICP-OES	BKK_EL0037	29-Feb-24	28-Feb-25	12
Water Lab	Magnesium	Hot Block	BKK_EL0054	22-Sep-23	22-Mar-25	18
Water Lab	Magnesium	Chamber (Cold Room)	BKK_EN0167	6-Dec-23	6-Jun-25	18
Water Lab	Copper	ICP-AES	BKK_EL0043	8-Oct-24	3-Apr-26	18
Water Lab	Copper	Hot Block	BKK_EL0054	22-Sep-23	22-Mar-25	18
Water Lab	Copper	Chamber (Cold Room)	BKK_EN0167	6-Dec-23	6-Jun-25	18
Water Lab	Cadmium	ICP-AES	BKK_EL0043	8-Oct-24	3-Apr-26	18
Water Lab	Cadmium	Hot Block	BKK_EL0054	22-Sep-23	22-Mar-25	18
Water Lab	Cadmium	Chamber (Cold Room)	BKK_EN0167	6-Dec-23	6-Jun-25	18

1

alsglobal.com



right solutions.  
right partner.

รายการเครื่องมือที่ใช้ในการวิเคราะห์ / ทดสอบ

Sample Name	Parameter	Equipment Name	ID No.	Calibrated Date	Next Cal	Freq. Calibrate (Months)
Water Lab	Zinc	ICP-AES	BKK_EL0043	8-Oct-24	3-Apr-26	18
Water Lab	Zinc	Hot Block	BKK_EL0054	22-Sep-23	22-Mar-25	18
Water Lab	Zinc	Chamber (Cold Room)	BKK_EN0167	6-Dec-23	6-Jun-25	18
Water Lab	Nickel	ICP-AES	BKK_EL0043	8-Oct-24	3-Apr-26	18
Water Lab	Nickel	Hot Block	BKK_EL0054	22-Sep-23	22-Mar-25	18
Water Lab	Nickel	Chamber (Cold Room)	BKK_EN0167	6-Dec-23	6-Jun-25	18
Water Lab	Selenium	ICP-AES	BKK_EL0043	8-Oct-24	3-Apr-26	18
Water Lab	Selenium	Hot Block	BKK_EL0054	22-Sep-23	22-Mar-25	18
Water Lab	Selenium	Chamber (Cold Room)	BKK_EN0167	6-Dec-23	6-Jun-25	18
Water Lab	Barium	ICP-AES	BKK_EL0043	8-Oct-24	3-Apr-26	18
Water Lab	Barium	Hot Block	BKK_EL0054	22-Sep-23	22-Mar-25	18
Water Lab	Barium	Chamber (Cold Room)	BKK_EN0167	6-Dec-23	6-Jun-25	18
Water Lab	Arsenic	ICP-AES	BKK_EL0043	8-Oct-24	3-Apr-26	18
Water Lab	Arsenic	Hot Block	BKK_EL0054	22-Sep-23	22-Mar-25	18
Water Lab	Arsenic	Chamber (Cold Room)	BKK_EN0167	6-Dec-23	6-Jun-25	18
Water Lab	Mercury	Mercury Analyzer	BKK_EL0128	6-Dec-23	6-Dec-25	12
Water Lab	Total Coliform	Autoclave	BKK_ML0037	17-Jul-23	17-Jan-25	18
Water Lab	Total Coliform	Incubator	BKK_ML0010	17-Jul-23	17-Jan-25	18
Water Lab	Total Coliform	Hot Air Oven	BKK_ML0013	23-Apr-24	23-Oct-25	18
Water Lab	Fecal Coliform	Autoclave	BKK_ML0037	17-Jul-23	17-Jan-25	18
Water Lab	Fecal Coliform	Incubator	BKK_ML0010	17-Jul-23	17-Jan-25	18
Water Lab	Fecal Coliform	Hot Air Oven	BKK_ML0013	23-Apr-24	23-Oct-25	18
Water Lab	Fecal Coliform	Water Bath	BKK_ML0056	1-Mar-24	1-Mar-25	12

2

alsglobal.com

## SITHIPORN ASSOCIATES CO., LTD. CALIBRATION LABORATORY

451-451/1 Sirthiporn Road, Bangbunmu, Bangkok, 10700 Thailand  
Tel. +66 2432 8331 Email: calibration@sithiporn.com

SITHIPORN  
ASSOCIATES



Cert. No. : ACC24009  
Pages : 1 of 3

## Calibration Certificate

Equipment : SOUND CALIBRATOR  
Manufacturer : RION  
Model : NC-74  
Serial No.: 34178118  
ID No.: BKK\_FS0631

Condition As Found : GOOD

Customer : ALS LABORATORY GROUP (THAILAND) CO., LTD.  
104 PHATTHANAKAN 40, PHATTHANAKAN ROAD,  
KHAENG PHATTHANAKAN, KHET SUAN LUANG,  
BANGKOK, 10250 THAILAND.

Location :  
Ambient Temperature : ( 23.0 ± 3 ) °C  
Pressure : ( 101.3 ± 3 ) kPa  
Relative Humidity : ( 50.0 ± 20 ) %

Received Date : 19 JANUARY 2024  
Calibration Date : 26 JANUARY 2024  
Date of Issue : 29 JANUARY 2024

Calibrated by : Nithakorn Pisutpaisan

Approved by :   
( Thanakul Petchurai )

This certificate is issued in accordance with the requirements of ISO/IEC 17025 standard, may not be reproduced other than in full, except with the prior written approval of the head of Calibration Laboratory.

## SITHIPORN ASSOCIATES CO., LTD. CALIBRATION LABORATORY

451-451/1 Sirthiporn Road, Bangbunmu, Bangkok, 10700 Thailand  
Tel. +66 2432 8331 Email: calibration@sithiporn.com

SITHIPORN  
ASSOCIATES



Cert. No. : ACC24009  
Job No. : VC67AC0059  
Pages : 2 of 3

Calibration Procedure : CP-AC-03

### Calibration Method :

This equipment was calibrated by follow on IEC-60942-2003 Standard.

The sound pressure level, frequency and total distortion of the sound calibrator was measured using the reference microphone.

### Condition of this result of calibration :

#### 1. Reference Standard Instruments :

Instrument	Model	Serial No.	Cert. No.	Due Date
Waveform Generator	33511B	MY52302742	EF-0010-23	07-FEB-24
Digital Multimeter	33461A	MY53220104	EELBP 30/0266	13-FEB-24
Digital Multimeter	33461A	MY53220076	EELBP 30/0267	13-FEB-24
Digital Multimeter	33461A	MY60024273	EELBP 31/0266	14-FEB-24
Programmable Attenuator	MAT-1070	62100114	EF-0011-23	08-FEB-24
Condenser Microphone	4180	2977900	AA-1001-23	14-FEB-24
Measuring Amplifier	NA-42KAI	34560495	AA-3002-23	14-FEB-24
Audio Analyzer	AVR-3360A	V744B6069	EF-0012-23	10-FEB-24

2. This result of calibration was found accurate as shown on date and place of calibration for this calibrated item only.

3. This certificate is traceable to the international system of unit maintained at :

3.1 National Institute of Metrology (Thailand).

3.2 Thailand Institute of Scientific and Technological Research (TISTR).

Cert. No. : ACC24009  
Job No. : VC67AC0059  
Pages : 3 of 3

## Result of calibration :

## 1. Sound pressure level

Specified sound pressure level (dB)	Measured value (dB)	Deviated value (dB)	Uncertainty (dB)	Acceptance limit (dB)
94	94.07	0.07	0.14	0.40

## 2. Frequency

Specified Frequency (Hz)	Measured value (Hz)	Deviated value (%)	Uncertainty (%)	Acceptance limit (%)
1000	1001.7	0.2	0.1	1.0

## 3. Total distortion

Measured value (%)	Uncertainty (%)	Acceptance limit (%)
1.65	0.10	3.0

The reported uncertainty is based on a standard uncertainty multiplied by coverage factor  $k = 2$  or any value following calculation, providing a level of confidence of approximately 95 %

End of Calibration Certificate

T. Petchuraj

Cert. No. : ACL23370  
Pages : 1 of 8

## Calibration Certificate

Equipment : SOUND LEVEL METER  
Manufacturer : RION  
Model : NL-42 / Microphone UC-52 / Preamplifier NH-24  
Serial No.: 01122547 / 143452 / 22584  
ID No.: BKK\_FS0034

Condition As Found : GOOD

Customer : ALS LABORATORY GROUP (THAILAND) CO., LTD.  
104 PHATTHANAKAN 40, PHATTHANAKAN ROAD,  
KHAENG PHATTHANAKAN, KHET SUAN LUANG,  
BANGKOK, 10250 THAILAND.

Location : -  
Ambient Temperature : ( 23.0 ± 3 ) °C  
Pressure : ( 101.3 ± 3 ) kPa  
Relative Humidity : ( 50.0 ± 20 ) %

Received Date : 07 NOVEMBER 2023  
Calibration Date : 29-30 NOVEMBER 2023  
Date of Issue : 06 DECEMBER 2023

Calibrated by : Nathakorn Pisutpaisan

Approved by :

T. Petchuraj  
( Thanakul Petchuraj )

This certificate is issued in accordance with the requirements of ISO/IEC 17025 standard, may not be reproduced other than in full, except with the prior written approval of the head of Calibration Laboratory.

QF-TS12-04-04-020664

## Continuation of Calibration Certificate

Cert. No. : ACL23370  
Job No. : VC67AC0025  
Pages : 2 of 8

Calibration Procedure : CP-AC-01

## Calibration Method :

This equipment was calibrated by based on IEC-61672-3 (2013) Standard for sound level meter (SLM).  
The SLM had tests to Acoustical and Electrical signal tests of frequency weighting with Anechoic chamber and Reference Standard Instruments.  
For tests results of each items were made by observation of each Instruments display and also with SLM's display.

## Condition of this result of calibration :

## 1. Reference Standard Instruments :

Instrument	Model	Serial No.	Cert. No.	Due Date
Waveform Generator	33210A	MY48017076	EF-0009-23	07-FEB-24
Waveform Generator	33511B	MY52302742	EF-0010-23	07-FEB-24
Digital Multimeter	33461A	MY53220104	EEL-BP 30/0266	13-FEB-24
Digital Multimeter	33461A	MY53220076	EEL-BP 29/0266	13-FEB-24
Digital Multimeter	34461A	MY60024273	EEL-BP 31/0266	14-FEB-24
Programmable Attenuator	MAT-1070	62100114	EF-0011-23	08-FEB-24
Condenser Microphone	4180	2977900	AA-1001-23	14-FEB-24
Measuring Amplifier	NA-42KAI	34560495	AA-3002-23	14-FEB-24

2. This result of calibration was found accurate as shown on date and place of calibration for this calibrated item only.

3. This certificate is traceable to the international system of unit maintained at :

3.1 National Institute of Metrology (Thailand).

3.2 Thailand Institute of Scientific and Technological Research (TISTR).

QF-TS12-04-04-020664

T. Petchuraj

## Continuation of Calibration Certificate

Cert. No. : ACL23370  
Job No. : VC67AC0025  
Pages : 3 of 8

## Summary of Measurement Result :

Parameter	Pass	Fail	Uncertainty (dB)	Maximum-permitted uncertainty of measurement (dB)
1. Absolute sensitivity	✓	-	0.2	N/A
2. Self-generated noise	✓	-	0.2	N/A
3. Acoustical signal tests of frequency weightings				
125 Hz	✓	-	0.3	0.6
1000 Hz	✓	-	0.3	0.6
8000 Hz	✓	-	0.3	0.7
4. Electrical signal tests of frequency weightings				
For 10 Hz to 4 kHz	✓	-	0.3	0.6
For > 4 kHz to 10 kHz	✓	-	0.3	0.7
For > 10 kHz to 20 kHz	-	-	-	1.0
5. Frequency and time weightings at 1 kHz	✓	-	0.2	0.2
6. Long - term stability	✓	-	0.1	0.1
7. Level linearity on the reference level range	✓	-	0.2	0.3
8. Level linearity including the level range control	✓	-	0.2	0.3
9. Tone burst response	✓	-	0.2	0.3
10. Peak C sound level	✓	-	0.2	0.35
11. Overload indication	✓	-	0.2	0.25
12. High level stability	✓	-	0.1	0.1

Note : Pass/Fail evaluation for each parameter, will be considered together from the acceptance limit and the Maximum-permitted uncertainty of measurement.

QF-TS12-04-04-020664

T. Petchuraj



## Continuation of Calibration Certificate

Cert. No. : ACL23370  
Job No. : VC67AC0025  
Pages : 4 of 8

## Result of calibration :

## 1. Absolute sensitivity

Reference Acoustic Signal ( dB )	Measured Value ( dB )	Deviation ( dB )	Acceptance Limit ( dB )
93.9 (93.98)	93.9	0.0	±0.3

## 2. Self-generated noise

## 2.1 Normal test

Measured Value ( dB )
15.6

## 2.2 The microphone of the sound level meter was replaced by electrical signal input device.

Frequency Weighting	Measured value ( dB )
A - weight	10.8
C - weight	17.0
Flat	22.9

## 3. Acoustical signal tests of frequency weightings

Meter free-field acoustic response at a level of 84 dB

Frequency ( Hz )	Deviation from various frequency weighting response curve (dB)			
	Flat	C-weight	A-weight	Acceptance Limits
125	0.9	1.0	0.9	± 1.5
1000	0.1	0.1	0.1	± 1.0
8000	-1.7	-1.7	-1.6	±5.0

QF-TS12-04-04-020664

T. Retin

## Continuation of Calibration Certificate

Cert. No. : ACL23370  
Job No. : VC67AC0025  
Pages : 5 of 8

## 4. Electrical signal tests of frequency weightings

Weighting network response with relative to 1 kHz.

Frequency ( Hz )	Deviation from various frequency weighting response curve (dB)			
	Flat	C-weight	A-weight	Acceptance Limits
63	0.0	-0.1	0.0	±2.0
125	0.0	0.0	0.0	±1.5
250	0.0	0.0	0.0	±1.5
500	0.0	0.0	0.0	±1.5
1000	0.0	0.0	0.0	±1.0
2000	0.0	0.0	0.0	±2.0
4000	0.0	0.0	0.0	±3.0
8000	0.0	0.1	0.1	±5.0

## 5. Frequency and time weightings at 1 kHz

## 5.1 Frequency weightings at 1 kHz

Frequency Weighting	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
A - weight	94.0	94.0	0.0	± 0.2
C - weight	94.0	94.0	0.0	± 0.2
Flat	94.0	94.0	0.0	± 0.2

## 5.2 Time weighting at 1 kHz

Frequency Weighting	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Fast	94.0	94.0	0.0	± 0.1
Slow	94.0	94.0	0.0	± 0.1
Leq	94.0	94.0	0.0	± 0.1

## 6. Long - term stability

Frequency Weighting	SLM Display at initial ( dB )	SLM Display at final ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
A - weight	94.0	94.0	0.0	± 0.3

QF-TS12-04-04-020664

T. Retin

## Continuation of Calibration Certificate

Cert. No. : ACL23370  
Job No. : VC67AC0025  
Pages : 6 of 8

## 7. Level linearity on the reference level range

Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
137.0	137.0	0.0	± 1.1
136.0	136.0	0.0	± 1.1
135.0	135.0	0.0	± 1.1
134.0	134.0	0.0	± 1.1
133.0	133.0	0.0	± 1.1
132.0	132.0	0.0	± 1.1
131.0	131.0	0.0	± 1.1
129.0	129.0	0.0	± 1.1
124.0	124.0	0.0	± 1.1
119.0	119.0	0.0	± 1.1
114.0	114.0	0.0	± 1.1
109.0	109.0	0.0	± 1.1
104.0	104.0	0.0	± 1.1
99.0	99.0	0.0	± 1.1
94.0	94.0	0.0	± 1.1
89.0	89.0	0.0	± 1.1
84.0	84.0	0.0	± 1.1
79.0	79.0	0.0	± 1.1
74.0	74.0	0.0	± 1.1
69.0	69.0	0.0	± 1.1
64.0	64.0	0.0	± 1.1
59.0	59.0	0.0	± 1.1
54.0	53.9	-0.1	± 1.1
49.0	49.0	0.0	± 1.1
44.0	44.0	0.0	± 1.1
39.0	38.9	-0.1	± 1.1
34.0	34.0	0.0	± 1.1
30.0	30.0	0.0	± 1.1
29.0	29.0	0.0	± 1.1
28.0	28.0	0.0	± 1.1
27.0	27.0	0.0	± 1.1
26.0	26.0	0.0	± 1.1
25.0	25.0	0.0	± 1.1

QF-TS12-04-04-020664

T. Retin

## Continuation of Calibration Certificate

Cert. No. : ACL23370  
Job No. : VC67AC0025  
Pages : 7 of 8

## 8. Level linearity including the level range control

Range	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Auto	94.0	94.0	0.0	±1.1

## 9. Tone burst response

Time Weighting	Tone burst duration, Tb ( ms )	Cycle	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Fast	0.25	1	108.0	107.9	-0.1	1.5 ; -5.0
	2	8	117.0	117.0	0.0	1.0 ; -2.5
	200	800	134.0	134.0	0.0	±1.0
Slow	2	8	108.0	108.0	0.0	1.5 ; -5.0
	200	800	127.6	127.6	0.0	±1.0
SEL	0.25	1	99.0	98.9	-0.1	1.5 ; -5.0
	2	8	108.0	108.0	0.0	1.0 ; -2.5
	200	800	128.0	128.0	0.0	±1.0

## 10. Peak C sound level

Number of cycle in test signal	Anticipated Value ( dB )	Measured Value, L <sub>peak</sub> ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Continuous	133.0	133.0	0.0	±3.0
One	136.4	136.2	-0.2	±3.0

Number of cycle in test signal	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Continuous	133.0	133.0	0.0	±2.0
Positive half cycle	135.4	135.2	-0.2	±2.0
Negative half cycle	135.4	135.2	-0.2	±2.0

QF-TS12-04-04-020664

T. Retin



## Continuation of Calibration Certificate

Cert. No. : ACL23370  
Job No. : VC67AC0025  
Pages : 8 of 8

## 11. Overload indication

Measured value ( dB )		Deviated Value ( dB )	Acceptance Limits ( dB )
Positive one-half cycle	Negative one-half cycle		
89.7	89.8	0.1	±1.5

## 12. High level stability

Frequency Weighting	SLM Display at initial ( dB )	SLM Display at final ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
A - weight	137.0	137.0	0.0	±0.3

The reported uncertainty is based on a standard uncertainty multiplied by coverage factor  $k = 2$   
or any value following calculation, providing a level of confidence of approximately 95 %

End of Calibration Certificate

QF-TS12-04-04-020664

T. Petchur

451-451/1 Sridinthorn Rd., Bangbunru, Bangplud Bangkok 10700 THAILAND  
Tel:0-2435-8800 Fax:0-2433-1679 e-mail:cal-center@sithiporn.com http://www.sithiporn.comCert. No. : ACL23265  
Pages : 1 of 8

## Calibration Certificate

Equipment : SOUND LEVEL METER  
Manufacturer : RION  
Model : NL-42/ Microphone UC-52 / Preamplifier NH-24  
Serial No. : 01022263 / 136951 / 22311  
ID No. : BKK\_FS0032

Condition As Found : GOOD

Customer : ALS LABORATORY GROUP (THAILAND) CO., LTD.  
104 PHATTHANAKAN 40, PHATTHANAKAN ROAD,  
KHWAENG PHATTHANAKAN, KHET SUAN LUANG,  
BANGKOK, 10250 THAILAND.

Location :  
Ambient Temperature : ( 23.0 ± 3 ) °C  
Pressure : ( 101.3 ± 3 ) kPa  
Relative Humidity : ( 50.0 ± 20 ) %

Received Date : 23 AUGUST 2023  
Calibration Date : 01 SEPTEMBER 2023  
Date of Issue : 04 SEPTEMBER 2023



Calibrated by : Nathakorn Pisutpaisan

Approved by :

T. Petchur  
( Thanakul Petchurai )

This certificate is issued in accordance with the requirements of ISO/IEC 17025 standard, may not be reproduced  
other than in full, except with the prior written approval of the head of Calibration Laboratory.

QF-TS12-04-04-020664

## Continuation of Calibration Certificate

Cert. No. : ACL23265  
Job No. : VC66AC0094  
Pages : 2 of 8

Calibration Procedure : CP-AC-01

## Calibration Method :

This equipment was calibrated by based on IEC-61672-3 (2013) Standard for sound level meter (SLM).  
The SLM had tests to Acoustical and Electrical signal tests of frequency weighting with Anechoic chamber and Reference  
Standard Instruments.

For tests results of each items were made by observation of each Instruments display and also with SLM's display.

## Condition of this result of calibration :

## 1. Reference Standard Instruments :

Instrument	Model	Serial No.	Cert. No.	Due Date
Waveform Generator	33210A	MY48017076	EF-0009-23	07-FEB-24
Waveform Generator	33511B	MY52302742	EF-0010-23	07-FEB-24
Digital Multimeter	33461A	MY53220104	EEL-BP 30/0266	13-FEB-24
Digital Multimeter	33461A	MY53220076	EEL-BP 29/0266	13-FEB-24
Digital Multimeter	34461A	MY60024273	EEL-BP 31/0266	14-FEB-24
Programmable Attenuator	MAT-1070	62100114	EF-0011-23	08-FEB-24
Condenser Microphone	4180	2977900	AA-1001-23	14-FEB-24
Measuring Amplifier	NA-42KA1	34560495	AA-3002-23	14-FEB-24

2. This result of calibration was found accurate as shown on date and place of calibration for this calibrated item only.

3. This certificate is traceable to the international system of unit maintained at :

3.1 National Institute of Metrology (Thailand).

3.2 Thailand Institute of Scientific and Technological Research (TISTR).

QF-TS12-04-04-020664

T. Petchur

## Continuation of Calibration Certificate

Cert. No. : ACL23265  
Job No. : VC66AC0094  
Pages : 3 of 8

## Summary of Measurement Result :

Parameter	Pass	Fail	Uncertainty (dB)	Maximum-permitted uncertainty of measurement (dB)
1. Absolute sensitivity	✓	-	0.2	N/A
2. Self-generated noise	✓	-	0.2	N/A
3. Acoustical signal tests of frequency weightings				
125 Hz	✓	-	0.3	0.6
1000 Hz	✓	-	0.3	0.6
8000 Hz	✓	-	0.3	0.7
4. Electrical signal tests of frequency weightings				
For 10 Hz to 4 kHz	✓	-	0.3	0.6
For > 4 kHz to 10 kHz	✓	-	0.3	0.7
For > 10 kHz to 20 kHz	-	-	-	1.0
5. Frequency and time weightings at 1 kHz	✓	-	0.2	0.2
6. Long - term stability	✓	-	0.1	0.1
7. Level linearity on the reference level range	✓	-	0.2	0.3
8. Level linearity including the level range control	✓	-	0.2	0.3
9. Tone burst response	✓	-	0.2	0.3
10. Peak C sound level	✓	-	0.2	0.35
11. Overload indication	✓	-	0.2	0.25
12. High level stability	✓	-	0.1	0.1

Note : Pass/Fail evaluation for each parameter,  
will be considered together from the acceptance limit and the Maximum-permitted uncertainty of measurement.

QF-TS12-04-04-020664

T. Petchur



## Continuation of Calibration Certificate

Cert. No. : ACL23265  
Job No. : VC66AC0094  
Pages : 4 of 8

## Result of calibration :

## 1. Absolute sensitivity

Reference Acoustic Signal ( dB )	Measured Value ( dB )	Deviation ( dB )	Acceptance Limit ( dB )
93.9 (93.98)	93.9	0.0	±0.3

## 2. Self-generated noise

## 2.1 Normal test

Measured Value ( dB )
16.7

## 2.2 The microphone of the sound level meter was replaced by electrical signal input device.

Frequency Weighting	Measured value ( dB )
A - weight	11.6
C - weight	17.6
Flat	23.4

## 3. Acoustical signal tests of frequency weightings

Meter free-field acoustic response at a level of 84 dB

Frequency ( Hz )	Deviation from various frequency weighting response curve (dB)			
	Flat	C-weight	A-weight	Acceptance Limits
125	0.5	0.5	0.5	± 1.5
1000	0.1	0.1	0.1	± 1.0
8000	-1.9	-1.8	-1.8	±5.0

QF-TS12-04-04-020664

7. Petch

## Continuation of Calibration Certificate

Cert. No. : ACL23265  
Job No. : VC66AC0094  
Pages : 5 of 8

## 4. Electrical signal tests of frequency weightings

Weighting network response with relative to 1 kHz.

Frequency ( Hz )	Deviation from various frequency weighting response curve (dB)			
	Flat	C-weight	A-weight	Acceptance Limits
63	-0.1	0.0	-0.1	±2.0
125	0.0	0.0	0.0	±1.5
250	0.0	0.0	0.0	±1.5
500	0.0	0.0	0.0	±1.5
1000	0.0	0.0	0.0	±1.0
2000	0.0	0.1	0.0	±2.0
4000	0.0	0.0	0.0	±3.0
8000	0.0	0.1	0.1	±5.0

## 5. Frequency and time weightings at 1 kHz

## 5.1 Frequency weightings at 1 kHz

Frequency Weighting	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
A - weight	94.0	94.0	0.0	± 0.2
C - weight	94.0	94.0	0.0	± 0.2
Flat	94.0	94.0	0.0	± 0.2

## 5.2 Time weighting at 1 kHz

Frequency Weighting	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Fast	94.0	94.0	0.0	± 0.1
Slow	94.0	94.0	0.0	± 0.1
Leq	94.0	94.0	0.0	± 0.1

## 6. Long - term stability

Frequency Weighting	SLM Display at initial ( dB )	SLM Display at final ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
A - weight	94.0	94.0	0.0	± 0.3

QF-TS12-04-04-020664

7. Petch

## Continuation of Calibration Certificate

Cert. No. : ACL23265  
Job No. : VC66AC0094  
Pages : 6 of 8

## 7. Level linearity on the reference level range

Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
137.0	137.0	0.0	± 1.1
136.0	136.0	0.0	± 1.1
135.0	135.0	0.0	± 1.1
134.0	134.0	0.0	± 1.1
133.0	133.0	0.0	± 1.1
132.0	132.0	0.0	± 1.1
131.0	131.0	0.0	± 1.1
129.0	129.0	0.0	± 1.1
124.0	124.0	0.0	± 1.1
119.0	119.0	0.0	± 1.1
114.0	114.0	0.0	± 1.1
109.0	109.0	0.0	± 1.1
104.0	104.0	0.0	± 1.1
99.0	99.0	0.0	± 1.1
94.0	94.0	0.0	± 1.1
89.0	89.0	0.0	± 1.1
84.0	84.0	0.0	± 1.1
79.0	79.0	0.0	± 1.1
74.0	74.0	0.0	± 1.1
69.0	69.0	0.0	± 1.1
64.0	64.0	0.0	± 1.1
59.0	59.0	0.0	± 1.1
54.0	54.0	0.0	± 1.1
49.0	49.0	0.0	± 1.1
44.0	44.0	0.0	± 1.1
39.0	39.0	0.0	± 1.1
34.0	34.0	0.0	± 1.1
30.0	29.9	-0.1	± 1.1
29.0	28.9	-0.1	± 1.1
28.0	27.9	-0.1	± 1.1
27.0	26.9	-0.1	± 1.1
26.0	25.9	-0.1	± 1.1
25.0	24.9	-0.1	± 1.1

QF-TS12-04-04-020664

7. Petch

## Continuation of Calibration Certificate

Cert. No. : ACL23265  
Job No. : VC66AC0094  
Pages : 7 of 8

## 8. Level linearity including the level range control

Range	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Auto	94.0	94.0	0.0	±1.1

## 9. Tone burst response

Time Weighting	Tone burst duration, Tb ( ms )	Cycle	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Fast	0.25	1	108.0	107.9	-0.1	1.5 ; -5.0
	2	8	117.0	117.0	0.0	1.0 ; -2.5
	200	800	134.0	134.1	0.1	±1.0
Slow	2	8	108.0	108.0	0.0	1.5 ; -5.0
	200	800	127.6	127.6	0.0	±1.0
SEL	0.25	1	99.0	98.9	-0.1	1.5 ; -5.0
	2	8	108.0	108.0	0.0	1.0 ; -2.5
	200	800	128.0	128.1	0.1	±1.0

## 10. Peak C sound level

Number of cycle in test signal	Anticipated Value ( dB )	Measured Value, L <sub>peak</sub> ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Continuous	133.0	133.0	0.0	±3.0
One	136.4	136.3	-0.1	±3.0

Number of cycle in test signal	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Continuous	133.0	133.0	0.0	±2.0
Positive half cycle	135.4	135.2	-0.2	±2.0
Negative half cycle	135.4	135.2	-0.2	±2.0

QF-TS12-04-04-020664

7. Petch



Continuation of Calibration Certificate

Cert. No. : ACL23265  
Job No. : VC66AC0094  
Pages : 8 of 8

11. Overload indication

Measured value ( dB )		Deviated Value ( dB )	Acceptance Limits ( dB )
Positive one-half cycle	Negative one-half cycle		
89.5	89.5	0.0	±1.5

12. High level stability

Frequency Weighting	SLM Display at initial ( dB )	SLM Display at final ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
A - weight	137.0	137.0	0.0	±0.3

The reported uncertainty is based on a standard uncertainty multiplied by coverage factor  $k = 2$  or any value following calculation, providing a level of confidence of approximately 95 %

End of Calibration Certificate

451-451/1 Srinthorn Road, Bangbunmu, Bangkok, 10700 Thailand  
Tel. +66 2433 8331 Email : calibration@sithiphorn.com

Cert. No. : ACL24226  
Pages : 1 of 8

Calibration Certificate

Equipment : SOUND LEVEL METER  
Manufacturer : RION  
Model : NL-42 / Microphone UC-52 / Preamplifier NH-24  
Serial No. : 00658241 / 158767 / 58769  
ID No. : BKK\_FS0098

Condition As Found : GOOD

Customer : ALS LABORATORY GROUP (THAILAND) CO., LTD.  
104 PHATTHANAKAN 40, PHATTHANAKAN ROAD,  
KHAENG PHATTHANAKAN, KHET SUAN LUANG,  
BANGKOK, 10250 THAILAND.

Location : -  
Ambient Temperature : ( 23.0 ± 3 ) °C  
Pressure : ( 101.3 ± 3 ) kPa  
Relative Humidity : ( 50.0 ± 20 ) %  
Received Date : 02 JULY 2024  
Calibration Date : 09-10 JULY 2024  
Date of Issue : 12 JULY 2024

REVIEW BY	Mahom P.
APPROVED BY	T. Petch.
NEXT CAL DATE	9/9/25

Calibrated by : Nathakorn Pisutpaisan

Approved by : T. Petch.  
( Thanakul Petchurai )

This certificate is issued in accordance with the requirements of ISO/IEC 17025 standard, may not be reproduced other than in full, except with the prior written approval of the head of Calibration Laboratory.

QF-TS12-04-04-020664

T. Petch.

451-451/1 Srinthorn Road, Bangbunmu, Bangkok, 10700 Thailand  
Tel. +66 2433 8331 Email : calibration@sithiphorn.com

Cert. No. : ACL24226  
Job No. : VC67AC0119  
Pages : 2 of 8

Calibration Procedure : CP-AC-01

Calibration Method :

This equipment was calibrated by follow on IEC-61672-3 (2013) Standard for sound level meter (SLM).  
The SLM had tests to Acoustical and Electrical signal tests of frequency weighting with Anechoic chamber and Reference Standard Instruments.  
For tests results of each items were made by observation of each Instruments display and also with SLM's display.

Condition of this result of calibration :

1. Reference Standard Instruments :

Instrument	Model	Serial No.	Cert. No.	Due Date
Waveform Generator	33210A	MY48017076	EF-0009-24	05-FEB-25
Waveform Generator	33511B	MY52302742	EF-0007-24	05-FEB-25
Digital Multimeter	33461A	MY53220104	EEL-BP 21/0267	13-FEB-25
Digital Multimeter	33461A	MY53220076	EEL-BP 20/0267	15-FEB-25
Digital Multimeter	34461A	MY60024273	EEL-BP 22/0267	15-FEB-25
Programmable Attenuator	MAT-1070	62100114	EF-0008-24	05-FEB-25
Condenser Microphone	4180	2977900	AA-1001-24	12-FEB-25
Measuring Amplifier	NA-42KAJ	34560495	AA-3001-24	05-FEB-25

2. This result of calibration was found accurate as shown on date and place of calibration for this calibrated item only.

3. This certificate is traceable to the international system of unit maintained at :

3.1 National Institute of Metrology (Thailand).

3.2 Thailand Institute of Scientific and Technological Research (TISTR).

T. Petch.

451-451/1 Srinthorn Road, Bangbunmu, Bangkok, 10700 Thailand  
Tel. +66 2433 8331 Email : calibration@sithiphorn.com

Cert. No. : ACL24226  
Job No. : VC67AC0119  
Pages : 3 of 8

Summary of Measurement Result :

Parameter	Uncertainty (dB)	Maximum-permitted uncertainty of measurement (dB)
1. Absolute sensitivity	0.2	N/A
2. Self-generated noise	0.2	N/A
3. Acoustical signal tests of frequency weightings		
125 Hz	0.3	0.6
1000 Hz	0.3	0.6
8000 Hz	0.3	0.7
4. Electrical signal tests of frequency weightings		
For 10 Hz to 4 kHz	0.3	0.6
For > 4 kHz to 10 kHz	0.3	0.7
For > 10 kHz to 20 kHz	-	1.0
5. Frequency and time weightings at 1 kHz	0.2	0.2
6. Long - term stability	0.1	0.1
7. Level linearity on the reference level range	0.2	0.3
8. Level linearity including the level range control	0.2	0.3
9. Tone burst response	0.2	0.3
10. Peak C sound level	0.2	0.35
11. Overload indication	0.2	0.25
12. High level stability	0.1	0.1

T. Petch.

Cert. No. : ACL24226  
Job No. : VC67AC0119  
Page : 4 of 8

**Result of calibration :**

**1. Absolute sensitivity**

Reference Acoustic Signal ( dB )	Measured Value ( dB )	Deviation ( dB )	Acceptance Limit ( dB )
93.9 (93.94)	93.9	0.0	±0.3

**2. Self-generated noise**

**2.1 Normal test**

Measured Value ( dB )
15.6

**2.2 The microphone of the sound level meter was replaced by electrical signal input device.**

Frequency Weighting	Weighting ( dB )
A - weight	11.2
C - weight	17.4
Flat	23.2

**3. Acoustical signal tests of frequency weightings**

Meter free-field acoustic response at a level of 84 dB

Frequency ( Hz )	Deviation from various frequency weighting response curve (dB)			
	Flat	C-weight	A-weight	Acceptance Limits
125	0.4	0.5	0.5	± 1.5
1000	-0.1	-0.1	-0.1	± 1.0
8000	-1.5	-1.4	-1.4	±5.0

*T. Petch.*

Cert. No. : ACL24226  
Job No. : VC67AC0119  
Pages : 5 of 8

**4. Electrical signal tests of frequency weightings**

Weighting network response with relative to 1 kHz.

Frequency ( Hz )	Deviation from various frequency weighting response curve (dB)			
	Flat	C-weight	A-weight	Acceptance Limits
63	-0.1	-0.1	-0.1	±2.0
125	0.0	0.0	-0.1	±1.5
250	0.0	0.0	0.0	±1.5
500	0.0	0.0	-0.1	±1.5
1000	0.0	0.0	0.0	±1.0
2000	0.0	0.0	0.0	±2.0
4000	0.0	0.0	0.0	±3.0
8000	0.0	0.1	0.1	±5.0

**5. Frequency and time weightings at 1 kHz**

**5.1 Frequency weightings at 1 kHz**

Frequency Weighting	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
A - weight	94.0	94.0	0.0	± 0.2
C - weight	94.0	94.0	0.0	± 0.2
Flat	94.0	94.0	0.0	± 0.2

**5.2 Time weighting at 1 kHz**

Frequency Weighting	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Fast	94.0	94.0	0.0	± 0.1
Slow	94.0	94.0	0.0	± 0.1
Leq	94.0	94.0	0.0	± 0.1

**6. Long - term stability**

Frequency Weighting	SLM Display at initial ( dB )	SLM Display at final ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
A - weight	94.0	94.0	0.0	± 0.3

*T. Petch.*

Cert. No. : ACL24226  
Job No. : VC67AC0119  
Pages : 6 of 8

**7. Level linearity on the reference level range**

Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
137.0	137.0	0.0	± 1.1
136.0	136.0	0.0	± 1.1
135.0	135.0	0.0	± 1.1
134.0	134.0	0.0	± 1.1
133.0	133.0	0.0	± 1.1
132.0	132.0	0.0	± 1.1
131.0	131.0	0.0	± 1.1
129.0	129.0	0.0	± 1.1
124.0	124.0	0.0	± 1.1
119.0	119.0	0.0	± 1.1
114.0	114.0	0.0	± 1.1
109.0	109.0	0.0	± 1.1
104.0	104.0	0.0	± 1.1
99.0	99.0	0.0	± 1.1
94.0	94.0	0.0	± 1.1
89.0	89.0	0.0	± 1.1
84.0	84.1	0.1	± 1.1
79.0	79.0	0.0	± 1.1
74.0	74.1	0.1	± 1.1
69.0	69.1	0.1	± 1.1
64.0	64.0	0.0	± 1.1
59.0	59.1	0.1	± 1.1
54.0	54.0	0.0	± 1.1
49.0	49.0	0.0	± 1.1
44.0	44.0	0.0	± 1.1
39.0	39.0	0.0	± 1.1
34.0	34.0	0.0	± 1.1
30.0	30.0	0.0	± 1.1
29.0	29.0	0.0	± 1.1
28.0	28.0	0.0	± 1.1
27.0	27.0	0.0	± 1.1
26.0	26.1	0.1	± 1.1
25.0	25.1	0.1	± 1.1

*T. Petch.*

Cert. No. : ACL24226  
Job No. : VC67AC0119  
Pages : 7 of 8

**8. Level linearity including the level range control**

Range	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Auto	94.0	94.0	0.0	±1.1

**9. Tone burst response**

Time Weighting	Tone burst duration, Tb ( ms )	Cycle	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Fast	0.25	1	108.0	107.9	-0.1	1.5 ; -5.0
	2	8	117.0	117.0	0.0	1.0 ; -2.5
	200	800	134.0	134.0	0.0	±1.0
Slow	2	8	108.0	108.0	0.0	1.5 ; -5.0
	200	800	127.6	127.6	0.0	±1.0
	0.25	1	99.0	98.9	-0.1	1.5 ; -5.0
SEL	2	8	108.0	108.0	0.0	1.0 ; -2.5
	200	800	128.0	128.0	0.0	±1.0

**10. Peak C sound level**

Number of cycle in test signal	Anticipated Value ( dB )	Measured Value, L <sub>peak</sub> ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Continuous	133.0	133.0	0.0	±3.0
One	136.4	136.4	0.0	±3.0

Number of cycle in test signal	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Continuous	133.0	133.1	0.1	±2.0
Positive half cycle	135.4	135.2	-0.2	±2.0
Negative half cycle	135.4	135.3	-0.1	±2.0

*T. Petch.*



# SITHIPORN ASSOCIATES CO., LTD. CALIBRATION LABORATORY

451-451/1 Srinthorn Road, Bangbunru, Bangkok, 10700 Thailand  
Tel. +66 2433 8331 Email: calibration@sithiporn.com



Cert. No. : ACI.24226  
Job No. : VC67AC0119  
Pages : 8 of 8

## 11. Overload indication

Measured value ( dB )		Deviated Value ( dB )	Acceptance Limits ( dB )
Positive one-half cycle	Negative one-half cycle		
89.5	89.7	0.2	±1.5

## 12. High level stability

Frequency Weighting	SLM Display at initial ( dB )	SLM Display at final ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
A - weight	137.0	137.0	0.0	±0.3

The reported uncertainty is based on a standard uncertainty multiplied by coverage factor  $k = 2$  or any value following calculation, providing a level of confidence of approximately 95 %

End of Calibration Certificate

*T. Petchur*

# SITHIPORN ASSOCIATES CO., LTD. CALIBRATION LABORATORY

451-451/1 Srinthorn Road, Bangbunru, Bangkok, 10700 Thailand  
Tel. +66 2433 8331 Email: calibration@sithiporn.com



Cert. No. : ACC24026  
Pages : 1 of 3

## Calibration Certificate

Equipment : SOUND CALIBRATOR  
Manufacturer : RION  
Model : NC-74  
Serial No. : 34178117  
ID No. : BKK\_FS0630

Condition As Found : GOOD

Customer : ALS LABORATORY GROUP (THAILAND) CO., LTD.  
104 PHATTHANAKAN 40, PHATTHANAKAN ROAD,  
KHWAENG PHATTHANAKAN, KHET SUAN LUANG,  
BANGKOK, 10250 THAILAND.

Location :  
Ambient Temperature : ( 23.0 ± 3 ) °C  
Pressure : ( 101.3 ± 3 ) kPa  
Relative Humidity : ( 50.0 ± 20 ) %

Received Date : 11 JUNE 2024  
Calibration Date : 27 JUNE 2024  
Date of Issue : 28 JUNE 2024



Calibrated by : Nathakorn Pisutpaisan

Approved by : *T. Petchur*  
( Thanakul Petchurai )

This certificate is issued in accordance with the requirements of ISO/IEC 17025 standard, may not be reproduced other than in full, except with the prior written approval of the head of Calibration Laboratory.

*T. Petchur*

# SITHIPORN ASSOCIATES CO., LTD. CALIBRATION LABORATORY

451-451/1 Srinthorn Road, Bangbunru, Bangkok, 10700 Thailand  
Tel. +66 2433 8331 Email: calibration@sithiporn.com



Cert. No. : ACC24026  
Job No. : VC67AC0105  
Pages : 2 of 3

Calibration Procedure : CP-AC-03

## Calibration Method :

This equipment was calibrated by follow on IEC-60942-2003 Standard.

The sound pressure level, frequency and total distortion of the sound calibrator was measured using the reference microphone.

## Condition of this result of calibration :

### 1. Reference Standard Instruments :

Instrument	Model	Serial No.	Cert. No.	Due Date
Waveform Generator	33511B	MY52302742	EP-0007-24	05-FEB-25
Digital Multimeter	33461A	MY53220104	EEL_BP 21/0267	13-FEB-25
Digital Multimeter	33461A	MY53220076	EEL_BP 20/0267	15-FEB-25
Digital Multimeter	33461A	MY60024273	EEL_BP 22/0267	15-FEB-25
Programmable Attenuator	MAT-1070	62100114	EF-0008-24	05-FEB-25
Condenser Microphone	4180	2977900	AA-1001-24	12-FEB-25
Measuring Amplifier	NA-42KA1	34560495	AA-3001-24	05-FEB-25
Audio Analyzer	AVR-3360A	V744B6069	EF-0009-24	09-FEB-25

2. This result of calibration was found accurate as shown on date and place of calibration for this calibrated item only.

3. This certificate is traceable to the international system of unit maintained at :

3.1 National Institute of Metrology (Thailand).

3.2 Thailand Institute of Scientific and Technological Research (TISTR).

*T. Petchur*

# SITHIPORN ASSOCIATES CO., LTD. CALIBRATION LABORATORY

451-451/1 Srinthorn Road, Bangbunru, Bangkok, 10700 Thailand  
Tel. +66 2433 8331 Email: calibration@sithiporn.com



Cert. No. : ACC24026  
Job No. : VC67AC0105  
Pages : 3 of 3

## Result of calibration :

### 1. Sound pressure level

Specified sound pressure level (dB)	Measured value (dB)	Deviated value (dB)	Uncertainty (dB)	Acceptance limit (dB)
94	94.14	0.14	0.14	0.40

### 2. Frequency

Specified Frequency (Hz)	Measured value (Hz)	Deviated value ( % )	Uncertainty ( % )	Acceptance limit ( % )
1000	1001.7	0.2	0.1	1.0

### 3. Total distortion

Measured value ( % )	Uncertainty ( % )	Acceptance limit ( % )
1.56	0.10	3.0

The reported uncertainty is based on a standard uncertainty multiplied by coverage factor  $k = 2$  or any value following calculation, providing a level of confidence of approximately 95 %

End of Calibration Certificate

*T. Petchur*

Cert. No. : ACL24262  
Pages : 1 of 8

## Calibration Certificate

Equipment : SOUND LEVEL METER  
Manufacturer : RION  
Model : NL-42 / Microphone UC-52 / Preamplifier NH-24  
Serial No.: 00658242 / 157782 / 48097  
ID No.: BKK\_FS0099

Condition As Found : GOOD

Customer : ALS LABORATORY GROUP (THAILAND) CO., LTD.  
104 PHATTHANAKAN 40, PHATTHANAKAN ROAD,  
KHWAENG PHATTHANAKAN, KHET SUAN LUANG,  
BANGKOK, 10250 THAILAND.

Location :  
Ambient Temperature : ( 23.0 ± 3 ) °C  
Pressure : ( 101.3 ± 3 ) kPa  
Relative Humidity : ( 50.0 ± 20 ) %

Received Date : 09 AUGUST 2024  
Calibration Date : 30 AUGUST 2024  
Date of Issue : 03 SEPTEMBER 2024



Calibrated by : Nathakorn Pisutpaisan

Approved by : *Thanakul P.*  
( Thanakul Petchurai )

This certificate is issued in accordance with the requirements of ISO/IEC 17025 standard, may not be reproduced  
other than in full, except with the prior written approval of the head of Calibration Laboratory.

Cert. No. : ACL24262  
Job No. : VC67AC0139  
Pages : 2 of 8

Calibration Procedure : CP-AC-01

## Calibration Method :

This equipment was calibrated by follow on IEC-61672-3 (2013) Standard for sound level meter (SLM).

The SLM had tests to Acoustical and Electrical signal tests of frequency weighting with Anechoic chamber and Reference Standard Instruments.

For tests results of each items were made by observation of each Instruments display and also with SLM's display.

## Condition of this result of calibration :

## 1. Reference Standard Instruments :

Instrument	Model	Serial No.	Cert. No.	Due Date
Waveform Generator	33210A	MY48017076	EF-0009-24	05-FEB-25
Waveform Generator	33511B	MY52302742	EF-0007-24	05-FEB-25
Digital Multimeter	33461A	MY53220104	EEL_BP 21/0267	13-FEB-25
Digital Multimeter	33461A	MY53220076	EEL_BP 20/0267	15-FEB-25
Digital Multimeter	34461A	MY60024273	EEL_BP 22/0267	15-FEB-25
Programmable Attenuator	MAT-1070	62100114	EF-0008-24	05-FEB-25
Condenser Microphone	4180	2977900	AA-1001-24	12-FEB-25
Measuring Amplifier	NA-42KAI	34560495	AA-3001-24	05-FEB-25

2. This result of calibration was found accurate as shown on date and place of calibration for this calibrated item only.

3. This certificate is traceable to the international system of unit maintained at :

- 3.1 National Institute of Metrology (Thailand),  
3.2 Thailand Institute of Scientific and Technological Research (TISTR).

*Thanakul P.*Cert. No. : ACL24262  
Job No. : VC67AC0139  
Pages : 3 of 8

## Summary of Measurement Result :

Parameter	Uncertainty (dB)	Maximum-permitted uncertainty of measurement (dB)
1. Absolute sensitivity	0.2	N/A
2. Self-generated noise	0.2	N/A
3. Acoustical signal tests of frequency weightings		
125 Hz	0.3	0.6
1000 Hz	0.3	0.6
8000 Hz	0.3	0.7
4. Electrical signal tests of frequency weightings		
For 10 Hz to 4 kHz	0.3	0.6
For > 4 kHz to 10 kHz	0.3	0.7
For > 10 kHz to 20 kHz	0.3	1.0
5. Frequency and time weightings at 1 kHz	0.2	0.2
6. Long - term stability	0.1	0.1
7. Level linearity on the reference level range	0.2	0.3
8. Level linearity including the level range control	0.2	0.3
9. Tone burst response	0.2	0.3
10. Peak C sound level	0.2	0.35
11. Overload indication	0.2	0.25
12. High level stability	0.1	0.1

*Thanakul P.*Cert. No. : ACL24262  
Job No. : VC67AC0139  
Page : 4 of 8

## Result of calibration :

## 1. Absolute sensitivity

Reference Acoustic Signal ( dB )	Measured Value ( dB )	Deviation ( dB )	Acceptance Limit ( dB )
93.9 (93.94)	93.9	0.0	±0.3

## 2. Self-generated noise

## 2.1 Normal test

Measured Value ( dB )
17.4

2.2 The microphone of the sound level meter was replaced by electrical signal input device.

Frequency Weighting ( dB )	Weighting ( dB )
A - weight	15.4
C - weight	21.4
Flat	27.2

## 3. Acoustical signal tests of frequency weightings

Meter free-field acoustic response at a level of 84 dB

Frequency ( Hz )	Deviation from various frequency weighting response curve (dB)			
	Flat	C-weight	A-weight	Acceptance Limits
125	0.1	0.2	0.2	± 1.5
1000	-0.1	-0.1	-0.1	± 1.0
8000	-1.5	-1.5	-1.4	±5.0

*Thanakul P.*



Cert. No. : ACL24262  
Job No. : VC67AC0139  
Pages : 5 of 8

#### 4. Electrical signal tests of frequency weightings

Weighting network response with relative to 1 kHz.

Frequency ( Hz )	Deviation from various frequency weighting response curve (dB)			
	Flat	C-weight	A-weight	Acceptance Limits
63	-0.1	0.0	0.0	±2.0
125	0.0	0.0	0.0	±1.5
250	0.0	0.0	0.0	±1.5
500	0.0	0.0	0.0	±1.5
1000	0.0	0.0	0.0	±1.0
2000	0.0	0.0	0.0	±2.0
4000	0.0	0.0	0.0	±3.0
8000	0.0	0.1	0.1	±5.0

#### 5. Frequency and time weightings at 1 kHz

##### 5.1 Frequency weightings at 1 kHz

Frequency Weighting	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
A - weight	94.0	94.0	0.0	± 0.2
C - weight	94.0	94.0	0.0	± 0.2
Flat	94.0	94.0	0.0	± 0.2

##### 5.2 Time weighting at 1 kHz

Frequency Weighting	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Fast	94.0	94.0	0.0	± 0.1
Slow	94.0	94.0	0.0	± 0.1
Leq	94.0	94.0	0.0	± 0.1

#### 6. Long - term stability

Frequency Weighting	SLM Display at initial ( dB )	SLM Display at final ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
A - weight	94.0	94.0	0.0	± 0.3

*T. Petch*

Cert. No. : ACL24262  
Job No. : VC67AC0139  
Pages : 6 of 8

#### 7. Level linearity on the reference level range

Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
137.0	137.0	0.0	± 1.1
136.0	136.0	0.0	± 1.1
135.0	135.0	0.0	± 1.1
134.0	134.0	0.0	± 1.1
133.0	133.0	0.0	± 1.1
132.0	132.0	0.0	± 1.1
131.0	131.0	0.0	± 1.1
129.0	129.0	0.0	± 1.1
124.0	124.0	0.0	± 1.1
119.0	119.0	0.0	± 1.1
114.0	114.0	0.0	± 1.1
109.0	109.0	0.0	± 1.1
104.0	104.0	0.0	± 1.1
99.0	99.0	0.0	± 1.1
94.0	94.0	0.0	± 1.1
89.0	89.1	0.1	± 1.1
84.0	84.1	0.1	± 1.1
79.0	79.0	0.0	± 1.1
74.0	74.1	0.1	± 1.1
69.0	69.1	0.1	± 1.1
64.0	64.0	0.0	± 1.1
59.0	59.1	0.1	± 1.1
54.0	54.0	0.0	± 1.1
49.0	49.0	0.0	± 1.1
44.0	44.0	0.0	± 1.1
39.0	39.0	0.0	± 1.1
34.0	34.1	0.1	± 1.1
30.0	30.2	0.2	± 1.1
29.0	29.3	0.3	± 1.1
28.0	28.4	0.4	± 1.1
27.0	27.5	0.5	± 1.1
26.0	26.6	0.6	± 1.1
25.0	25.7	0.7	± 1.1

*T. Petch*

Cert. No. : ACL24262  
Job No. : VC67AC0139  
Pages : 7 of 8

#### 8. Level linearity including the level range control

Range	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
130	94.0	94.0	0.0	±1.1

Range	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
130	30.0	30.1	0.1	±1.1

#### 9. Tone burst response

Time Weighting	Tone burst duration, Tb ( ms )	Cycle	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Fast	0.25	1	108.0	107.9	-0.1	1.5 ; -5.0
	2	8	117.0	117.0	0.0	1.0 ; -2.5
	200	800	134.0	134.0	0.0	±1.0
Slow	2	8	108.0	108.0	0.0	1.5 ; -5.0
	200	800	127.6	127.6	0.0	±1.0
SEL	0.25	1	99.0	98.9	-0.1	1.5 ; -5.0
	2	8	108.0	108.0	0.0	1.0 ; -2.5
	200	800	128.0	128.0	0.0	±1.0

*T. Petch*

Cert. No. : ACL24262  
Job No. : VC67AC0139  
Pages : 8 of 8

#### 10. Peak C sound level

Number of cycle in test signal	Anticipated Value ( dB )	Measured Value, L <sub>peak</sub> ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Continuous	133.0	133.0	0.0	±3.0
One	136.4	136.3	-0.1	±3.0

Number of cycle in test signal	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Continuous	133.0	133.0	0.0	±2.0
Positive half cycle	135.4	135.2	-0.2	±2.0
Negative half cycle	135.4	135.2	-0.2	±2.0

#### 11. Overload indication

Measured value ( dB )		Deviated Value ( dB )	Acceptance Limits ( dB )
Positive one-half cycle	Negative one-half cycle		
89.8	89.5	-0.3	±1.5

#### 12. High level stability

Frequency Weighting	SLM Display at initial ( dB )	SLM Display at final ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
A - weight	137.0	137.0	0.0	±0.3

The reported uncertainty is based on a standard uncertainty multiplied by coverage factor  $k = 2$   
or any value following calculation, providing a level of confidence of approximately 95 %

End of Calibration Certificate

*T. Petch*



Cert. No. : ACL24263  
Pages : 1 of 8

## Calibration Certificate


Equipment : SOUND LEVEL METER  
Manufacturer : RION  
Model : NL-42 / Microphone UC-52 / Preamplifier NH-24  
Serial No. : 00658243 / 157783 / 48098  
ID No. : BKK\_FS0100

Condition As Found : GOOD

Customer : ALS LABORATORY GROUP (THAILAND) CO., LTD.  
104 PHATTHANAKAN 40, PHATTHANAKAN ROAD,  
KHWAENG PHATTHANAKAN, KHET SUAN LUANG,  
BANGKOK, 10250 THAILAND.

Location : -  
Ambient Temperature : ( 23.0 ± 3 ) °C  
Pressure : ( 101.3 ± 3 ) kPa  
Relative Humidity : ( 50.0 ± 20 ) %

Received Date : 09 AUGUST 2024  
Calibration Date : 30 AUGUST 2024  
Date of Issue : 03 SEPTEMBER 2024

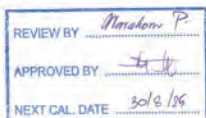

Calibrated by : Nathakorn Pisutpaisan

Approved by :

*T. Petchur*  
( Thanakul Petchurai )

This certificate is issued in accordance with the requirements of ISO/IEC 17025 standard, may not be reproduced  
other than in full, except with the prior written approval of the head of Calibration Laboratory.

Cert. No. : ACL24263  
Job No. : VC67AC0139  
Pages : 2 of 8

Calibration Procedure : CP-AC-01

## Calibration Method :

This equipment was calibrated by follow on IEC-61672-3 (2013) Standard for sound level meter (SLM).  
The SLM had tests to Acoustical and Electrical signal tests of frequency weighting with Anechoic chamber and Reference  
Standard Instruments.

For tests results of each items were made by observation of each Instruments display and also with SLM's display.

## Condition of this result of calibration :

## 1. Reference Standard Instruments :

Instrument	Model	Serial No.	Cert. No.	Due Date
Waveform Generator	33210A	MY48017076	EF-0009-24	05-FEB-25
Waveform Generator	33511B	MY52302742	EF-0007-24	05-FEB-25
Digital Multimeter	33461A	MY53220104	EELBP 21/0267	13-FEB-25
Digital Multimeter	33461A	MY53220076	EELBP 20/0267	15-FEB-25
Digital Multimeter	34461A	MY60024273	EELBP 22/0267	15-FEB-25
Programmable Attenuator	MAT-1070	62100114	EF-0008-24	05-FEB-25
Condenser Microphone	4180	2977900	AA-1001-24	12-FEB-25
Measuring Amplifier	NA-42KAI	34560495	AA-3001-24	05-FEB-25

2. This result of calibration was found accurate as shown on date and place of calibration for this calibrated item only.

3. This certificate is traceable to the international system of unit maintained at :

3.1 National Institute of Metrology (Thailand).

3.2 Thailand Institute of Scientific and Technological Research (TISTR).

*T. Petchur*Cert. No. : ACL24263  
Job No. : VC67AC0139  
Pages : 3 of 8

## Summary of Measurement Result :

Parameter	Uncertainty (dB)	Maximum-permitted uncertainty of measurement (dB)
1. Absolute sensitivity	0.2	N/A
2. Self-generated noise	0.2	N/A
3. Acoustical signal tests of frequency weightings		
125 Hz	0.3	0.6
1000 Hz	0.3	0.6
8000 Hz	0.3	0.7
4. Electrical signal tests of frequency weightings		
For 10 Hz to 4 kHz	0.3	0.6
For > 4 kHz to 10 kHz	0.3	0.7
For > 10 kHz to 20 kHz	0.3	1.0
5. Frequency and time weightings at 1 kHz	0.2	0.2
6. Long - term stability	0.1	0.1
7. Level linearity on the reference level range	0.2	0.3
8. Level linearity including the level range control	0.2	0.3
9. Tone burst response	0.2	0.3
10. Peak C sound level	0.2	0.35
11. Overload indication	0.2	0.25
12. High level stability	-0.1	0.1

*T. Petchur*Cert. No. : ACL24263  
Job No. : VC67AC0139  
Page : 4 of 8

## Result of calibration :

## 1. Absolute sensitivity

Reference Acoustic Signal ( dB )	Measured Value ( dB )	Deviation ( dB )	Acceptance Limit ( dB )
93.9 (93.94)	93.9	0.0	±0.3

## 2. Self-generated noise

## 2.1 Normal test

Measured Value ( dB )
15.4

2.2 The microphone of the sound level meter was replaced by electrical signal input device.

Frequency Weighting	Weighting ( dB )
A - weight	12.6
C - weight	18.8
Flat	24.4

## 3. Acoustical signal tests of frequency weightings

Meter free-field acoustic response at a level of 84 dB

Frequency ( Hz )	Deviation from various frequency weighting response curve (dB)			
	Flat	C-weight	A-weight	Acceptance Limits
125	0.2	0.2	0.2	±1.5
1000	0.0	0.0	0.0	±1.0
8000	-1.8	-1.7	-1.7	±5.0

*T. Petchur*



Cert. No. : ACL24263  
Job No. : VC67AC0139  
Pages : 5 of 8

## 4. Electrical signal tests of frequency weightings

Weighting network response with relative to 1 kHz.

Frequency ( Hz )	Deviation from various frequency weighting response curve (dB)			
	Flat	C-weight	A-weight	Acceptance Limits
63	0.0	-0.1	-0.1	±2.0
125	0.0	0.0	0.0	±1.5
250	0.0	0.0	-0.1	±1.5
500	0.0	0.0	-0.1	±1.5
1000	0.0	0.0	0.0	±1.0
2000	0.0	0.0	0.0	±2.0
4000	0.0	0.0	0.0	±3.0
8000	0.0	0.1	0.1	±5.0

## 5. Frequency and time weightings at 1 kHz

## 5.1 Frequency weightings at 1 kHz

Frequency Weighting	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
A - weight	94.0	94.0	0.0	± 0.2
C - weight	94.0	94.0	0.0	± 0.2
Flat	94.0	94.0	0.0	± 0.2

## 5.2 Time weighting at 1 kHz

Frequency Weighting	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Fast	94.0	94.0	0.0	± 0.1
Slow	94.0	94.0	0.0	± 0.1
Leq	94.0	94.0	0.0	± 0.1

## 6. Long - term stability

Frequency Weighting	SLM Display at initial ( dB )	SLM Display at final ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
A - weight	94.0	94.0	0.0	± 0.3

T. Petch

Cert. No. : ACL24263  
Job No. : VC67AC0139  
Pages : 6 of 8

## 7. Level linearity on the reference level range

Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
137.0	137.0	0.0	± 1.1
136.0	136.0	0.0	± 1.1
135.0	135.0	0.0	± 1.1
134.0	134.0	0.0	± 1.1
133.0	132.9	-0.1	± 1.1
132.0	131.9	-0.1	± 1.1
131.0	130.9	-0.1	± 1.1
129.0	129.0	0.0	± 1.1
124.0	124.0	0.0	± 1.1
119.0	119.0	0.0	± 1.1
114.0	114.0	0.0	± 1.1
109.0	109.0	0.0	± 1.1
104.0	104.0	0.0	± 1.1
99.0	99.0	0.0	± 1.1
94.0	94.0	0.0	± 1.1
89.0	89.0	0.0	± 1.1
84.0	84.0	0.0	± 1.1
79.0	79.0	0.0	± 1.1
74.0	74.0	0.0	± 1.1
69.0	69.0	0.0	± 1.1
64.0	64.0	0.0	± 1.1
59.0	59.1	0.1	± 1.1
54.0	54.0	0.0	± 1.1
49.0	49.0	0.0	± 1.1
44.0	44.0	0.0	± 1.1
39.0	39.1	0.1	± 1.1
34.0	34.2	0.2	± 1.1
30.0	30.6	0.6	± 1.1
29.0	29.6	0.6	± 1.1
28.0	28.7	0.7	± 1.1
27.0	27.9	0.9	± 1.1
26.0	27.3	1.3	± 1.1
25.0	26.4	1.4	± 1.1

T. Petch

Cert. No. : ACL24263  
Job No. : VC67AC0139  
Pages : 7 of 8

## 8. Level linearity including the level range control

Range	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
130	94.0	94.0	0.0	±1.1

Range	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
130	30.0	30.3	0.3	±1.1

## 9. Tone burst response

Time Weighting	Tone burst duration, Tb ( ms )	Cycle	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Fast	0.25	1	108.0	107.9	-0.1	1.5 ; -5.0
	2	8	117.0	117.0	0.0	1.0 ; -2.5
	200	800	134.0	134.0	0.0	±1.0
Slow	2	8	108.0	108.0	0.0	1.5 ; -5.0
	200	800	127.6	127.6	0.0	±1.0
SEL	0.25	1	99.0	98.9	-0.1	1.5 ; -5.0
	2	8	108.0	108.0	0.0	1.0 ; -2.5
	200	800	128.0	128.0	0.0	±1.0

T. Petch

Cert. No. : ACL24263  
Job No. : VC67AC0139  
Pages : 8 of 8

## 10. Peak C sound level

Number of cycle in test signal	Anticipated Value ( dB )	Measured Value, L <sub>peak</sub> ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Continuous	133.0	133.0	0.0	±3.0
One	136.4	135.5	-0.9	±3.0

Number of cycle in test signal	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Continuous	133.0	133.0	0.0	±2.0
Positive half cycle	135.4	135.1	-0.3	±2.0
Negative half cycle	135.4	135.1	-0.3	±2.0

## 11. Overload indication

Measured value ( dB )		Deviated Value ( dB )	Acceptance Limits ( dB )
Positive one-half cycle	Negative one-half cycle		
89.6	89.6	0.0	±1.5

## 12. High level stability

Frequency Weighting	SLM Display at initial ( dB )	SLM Display at final ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
A - weight	137.0	137.0	0.0	±0.3

The reported uncertainty is based on a standard uncertainty multiplied by coverage factor  $k = 2$   
or any value following calculation, providing a level of confidence of approximately 95 %

End of Calibration Certificate

T. Petch



Cert. No. : ACL24151  
Pages : 1 of 8

## Calibration Certificate

Equipment : SOUND LEVEL METER  
Manufacturer : RION  
Model : NL-42 / Microphone UC-52 / Preamplifier NH-24  
Serial No. : 00658244 / 158766 / 58768  
ID No. : BKK\_FS0101

Condition As Found : GOOD

Customer : ALS LABORATORY GROUP (THAILAND) CO., LTD.  
104 PHATTHANAKAN 40, PHATTHANAKAN ROAD,  
KHWAENG PHATTHANAKAN, KHET SUAN LUANG,  
BANGKOK, 10250 THAILAND.

Location : -  
Ambient Temperature : ( 23.0 ± 3 ) °C  
Pressure : ( 101.3 ± 3 ) kPa  
Relative Humidity : ( 50.0 ± 20 ) %

Received Date : 14 MAY 2024  
Calibration Date : 29 MAY 2024  
Date of Issue : 30 MAY 2024

Calibrated by : Nathakorn Pisurpaisan

Approved by : *[Signature]*  
( Thanakul Peichurai )

This certificate is issued in accordance with the requirements of ISO/IEC 17025 standard, may not be reproduced other than in full, except with the prior written approval of the head of Calibration Laboratory.

Cert. No. : ACL24151  
Job No. : VC67AC0089  
Pages : 2 of 8

Calibration Procedure : CP-AC-01

## Calibration Method :

This equipment was calibrated by follow on IEC-61672-3 (2013) Standard for sound level meter (SLM).  
The SLM had tests to Acoustical and Electrical signal tests of frequency weighting with Anechoic chamber and Reference Standard Instruments.

For tests results of each items were made by observation of each Instruments display and also with SLM's display.

## Condition of this result of calibration :

## 1. Reference Standard Instruments :

Instrument	Model	Serial No.	Cert. No.	Due Date
Waveform Generator	33210A	MY48017076	EF-0009-4	05-FEB-25
Waveform Generator	33511B	MY52302742	EF-0007-24	05-FEB-25
Digital Multimeter	33461A	MY53220104	EEL-BP 21/0267	13-FEB-25
Digital Multimeter	33461A	MY53220076	EEL-BP 20/0267	15-FEB-25
Digital Multimeter	34461A	MY60024273	EEL-BP 22/0267	15-FEB-25
Programmable Attenuator	MAT-1070	62100114	EF-0008-24	05-FEB-25
Condenser Microphone	4180	2977900	AA-1001-24	12-FEB-25
Measuring Amplifier	NA-42KAI	34560495	AA-3001-24	05-FEB-25

2. This result of calibration was found accurate as shown on date and place of calibration for this calibrated item only.

3. This certificate is traceable to the international system of unit maintained at :

- 3.1 National Institute of Metrology (Thailand).  
3.2 Thailand Institute of Scientific and Technological Research (TISTR).

Cert. No. : ACL24151  
Job No. : VC67AC0089  
Pages : 3 of 8

## Summary of Measurement Result :

Parameter	Uncertainty (dB)	Maximum-permitted uncertainty of measurement (dB)
1. Absolute sensitivity	0.2	N/A
2. Self-generated noise	0.2	N/A
3. Acoustical signal tests of frequency weightings		
125 Hz	0.3	0.6
1000 Hz	0.3	0.6
8000 Hz	0.3	0.7
4. Electrical signal tests of frequency weightings		
For 10 Hz to 4 kHz	0.3	0.6
For > 4 kHz to 10 kHz	0.3	0.7
For > 10 kHz to 20 kHz	-	1.0
5. Frequency and time weightings at 1 kHz	0.2	0.2
6. Long-term stability	0.1	0.1
7. Level linearity on the reference level range	0.2	0.3
8. Level linearity including the level range control	0.2	0.3
9. Tone burst response	0.2	0.3
10. Peak C sound level	0.2	0.35
11. Overload indication	0.2	0.25
12. High level stability	0.1	0.1

Cert. No. : ACL24151  
Job No. : VC67AC0089  
Pages : 4 of 8

## Result of calibration :

## 1. Absolute sensitivity

Reference Acoustic Signal (dB)	Measured Value (dB)	Deviation (dB)	Acceptance Limit (dB)
93.9 (93.94)	93.9	0.0	±0.3

## 2. Self-generated noise

## 2.1 Normal test

Measured Value (dB)
17.4

2.2 The microphone of the sound level meter was replaced by electrical signal input device.

Frequency Weighting	Measured value (dB)
A-weight	13.8
C-weight	20.5
Flat	25.9

## 3. Acoustical signal tests of frequency weightings

Meter free-field acoustic response at a level of 84 dB

Frequency (Hz)	Deviation from various frequency weighting response curve (dB)			
	Flat	C-weight	A-weight	Acceptance Limits
125	0.2	0.3	0.3	± 1.5
1000	0.0	0.0	0.0	± 1.0
8000	-0.5	-0.4	-0.4	± 5.0



4. Electrical signal tests of frequency weightings

Weighting network response with relative to 1 kHz.

Frequency ( Hz )	Deviation from various frequency weighting response curve (dB)			
	Flat	C-weight	A-weight	Acceptance Limits
63	-0.1	-0.1	0.0	±2.0
125	-0.1	0.0	0.0	±1.5
250	0.0	0.0	-0.1	±1.5
500	0.0	0.0	-0.1	±1.5
1000	0.0	0.0	0.0	±1.0
2000	0.0	0.0	0.0	±2.0
4000	0.0	0.0	0.0	±3.0
8000	0.0	0.1	0.1	±5.0

5. Frequency and time weightings at 1 kHz

5.1 Frequency weightings at 1 kHz

Frequency Weighting	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
A - weight	94.0	94.0	0.0	± 0.2
C- weight	94.0	94.0	0.0	± 0.2
Flat	94.0	94.0	0.0	± 0.2

5.2 Time weighting at 1 kHz

Frequency Weighting	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Fast	94.0	94.0	0.0	± 0.1
Slow	94.0	94.0	0.0	± 0.1
Leq	94.0	94.0	0.0	± 0.1

6. Long - term stability

Frequency Weighting	SLM Display at initial ( dB )	SLM Display at final ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
A - weight	94.0	94.0	0.0	± 0.3

*T. Petch*

7. Level linearity on the reference level range

Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
137.0	137.0	0.0	± 1.1
136.0	136.0	0.0	± 1.1
135.0	135.0	0.0	± 1.1
134.0	134.0	0.0	± 1.1
133.0	133.0	0.0	± 1.1
132.0	132.0	0.0	± 1.1
131.0	131.0	0.0	± 1.1
129.0	129.0	0.0	± 1.1
124.0	124.0	0.0	± 1.1
119.0	119.0	0.0	± 1.1
114.0	114.0	0.0	± 1.1
109.0	109.0	0.0	± 1.1
104.0	104.0	0.0	± 1.1
99.0	99.0	0.0	± 1.1
94.0	94.0	0.0	± 1.1
89.0	89.1	0.1	± 1.1
84.0	84.1	0.1	± 1.1
79.0	79.0	0.0	± 1.1
74.0	74.1	0.1	± 1.1
69.0	69.0	0.0	± 1.1
64.0	64.0	0.0	± 1.1
59.0	59.1	0.1	± 1.1
54.0	54.1	0.1	± 1.1
49.0	49.0	0.0	± 1.1
44.0	44.0	0.0	± 1.1
39.0	39.0	0.0	± 1.1
34.0	34.1	0.1	± 1.1
30.0	30.1	0.1	± 1.1
29.0	29.2	0.2	± 1.1
28.0	28.3	0.3	± 1.1
27.0	27.4	0.4	± 1.1
26.0	26.5	0.5	± 1.1
25.0	25.5	0.5	± 1.1

*T. Petch*

8. Level linearity including the level range control

Range	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Auto	94.0	94.0	0.0	±1.1

9. Tone burst response

Time Weighting	Tone burst duration, Tb ( ms )	Cycle	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Fast	0.25	1	108.0	107.9	-0.1	1.5 ; -5.0
	2	8	117.0	117.0	0.0	1.0 ; -2.5
	200	800	134.0	134.0	0.0	±1.0
Slow	2	8	108.0	108.0	0.0	1.5 ; -5.0
	200	800	127.6	127.6	0.0	±1.0
	200	800	128.0	128.1	0.1	±1.0

10. Peak C sound level

Number of cycle in test signal	Anticipated Value ( dB )	Measured Value, Lcpeak ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Continuous	133.0	133.0	0.0	±3.0
One	136.4	136.2	-0.2	±3.0

Number of cycle in test signal	Anticipated Value ( dB )	Measured Value ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
Continuous	133.0	133.0	0.0	±2.0
Positive half cycle	135.4	135.1	-0.3	±2.0
Negative half cycle	135.4	135.1	-0.3	±2.0

*T. Petch*

11. Overload indication

Measured value ( dB )		Deviated Value ( dB )	Acceptance Limits ( dB )
Positive one-half cycle	Negative one-half cycle		
89.6	89.6	0.0	±1.5

12. High level stability

Frequency Weighting	SLM Display at initial ( dB )	SLM Display at final ( dB )	Deviated Value ( dB )	Acceptance Limits ( dB )
A - weight	137.0	137.0	0.0	±0.3

The reported uncertainty is based on a standard uncertainty multiplied by coverage factor  $k = 2$   
or any value following calculation, providing a level of confidence of approximately 95 %

End of Calibration Certificate

*T. Petch*



## Certificate of Calibration

Cert.No.: 24CH1295  
Page.: 1 of 3

Equipment : pH Meter  
Manufacturer : Hach  
Model : HQ411d  
Serial No. : 200100031163  
ID No. : BKK\_EN0342  
Condition As-Received: Used Item  
Received Date : 16 October 2024  
Calibration Date : 17 October 2024  
Reference : 2410-0548DSC-5  
Submitted by : ALS Laboratory Group (Thailand) Co.,Ltd.  
104 Phatthanakan 40, Phatthanakan Rd.,  
Khwaeng Phatthanakan, Khet Suan Luang,  
Bangkok 10250 Thailand

Ambient Temperature : (25 ± 2.5) °C  
Relative Humidity : (50 ± 15) %  
Calibration Procedure : In - house method :  
- CP-CH5 by direct measurement with  
certified reference material (CRM)  
- CP-CH8 by comparison with temperature standard

Calibrated by : Warakorn Lerngatrakul

Approved by :

- ( ) Unnopphol Harachai  
( ) Ponpan Paipim  
(✓) Salthip Meangmai

Issue Date : 21 October 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



## Certificate of Calibration

Cert.No.: 24CG952  
Page.: 1 of 2

Equipment : Burette  
Capacity : 50 mL  
Serial No. : -  
ID. No. : BKK\_EN0171  
Manufacturer : Witeg  
Made in : Germany  
Submitted by : ALS Laboratory Group (Thailand) Co.,Ltd.  
104 Phatthanakan 40, Phatthanakan Rd.,  
Khwaeng Phatthanakan, Khet Suan Luang,  
Bangkok 10250 Thailand

Ambient Temperature : (20 ± 2.5) °C  
Relative Humidity : (50 ± 10) %  
Barometric Pressure : 760 mmHg  
Calibration Procedure : ASTM E 542 - 01

Calibrated by : Natcha Chayingcheiw

Approved by :

- ( ) Unnopphol Harachai  
(✓) Srisuda Khamtha  
( ) Sa-ngeunkam Wongsa

Issue Date : 27 February 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

### Calibration Results

Function : Temperature Measurement

(\*) Without adjustment

This equipment was connected with Temperature Probe;

- Model : PHC281

- Serial No. : 230473042902

Dimension of probe

- Length : 103 mm.

- Diameter : 12 mm.

- Immersion Depth : 90 mm.

Calibration Point (°C)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty of measurement (± °C)	Coverage factor k
25.0	25.002	25.0	-0.002	0.13	2.00

Remark : UUC\* = Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-oOo-





Equipment : Burette  
Received Date : 23 February 2024  
Condition As-Received : New Item  
Calibration Date : 27 February 2024  
Reference : 2402-0757DSC-1

Cert.No.: 24CG952  
Page.: 2 of 2

#### Condition of this result of calibration

##### 1. Reference Standard Instruments :

Instruments	Model	Serial No.	ID. No.	Certificate No.	Traceability	Due date
1) Balance	XP205DR	1126143764	140RC004	23MM538	TPA	15 Sep 2024
2) Thermo-Hygrograph	THDX-CE	00016540	140EC001	23H1275	TPA	09 June 2024
3) Thermometer	-	0834181	140EC005	23I948	TPA	10 Aug 2024

This certification is traceable to SI Unit

2. The certificate is valid only to the item calibrated on date and place of calibration.

3. True value is converted to true volume at the standard temperature of 20 °C

#### Calibration result :

Nominal capacity ( mL )	Reading ( mL )	Uncertainty ( ± mL )	k Factor
50	50.0032	0.010	2.00

**Remark** mL = cm<sup>3</sup>

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-o-o-



## Metrology

SCI ECO Services Company Limited

33/2 Moo 3, T.Banpa, A.Kaengkhoh, Saraburi 18110, Thailand.  
Saraburi Tel : +66 3627 3096 Fax : +66 3627 3100  
Bangkok Tel : +668 9205 6851 , +669 8247 2360  
Website : www.scieco.co.th E-Mail : calibrate@scg.com



Certificate No. T232160

Page 1 of 4

### Certificate of Calibration

Equipment : Chamber ( Cooling Room )

Manufacturer : KOLDTECH

Model : KM 320

Serial No. : TBN-1012061/05

Customer Code : BKK\_EN0167

ID No. : T2463A3

Customer : ALS Laboratory Group (Thailand) Co.,Ltd.


104 Phatthanakan 40, Phatthanakan Rd., Khwaeng Phatthanakan,

Khet Suan Luang, Bangkok 10250

Customer Location : Laboratory

Date of Receipt : 29 November 2023

Calibrated By : Atiphong Rongrat ( Technician )

Approved By :  / Boonchai Suriyawong (Site Calibration Manager)

Date of Issue : 09 JAN 2024

The uncertainties are for a confidence probability of approximately 95%.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standard laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the Metrology.

FM-L15 118/18-08-66



## Metrology

SCI ECO Services Company Limited

33/2 Moo 3, T.Banpa, A.Kaengkhoh, Saraburi 18110, Thailand.



Certificate No. T232160

Page 2 of 4

### Calibration Report

Equipment : Chamber ( Cooling Room )  
Date of Calibration : 6 December 2023  
Environment : Temperature : 23.4-24.9 °C  
Line Voltage : 221.4-230.2 V  
Relative Humidity : 55 - 65 %RH

#### Condition of this results of calibration :

1. This equipment was calibrated by insert 16 standard thermocouples type T into its chamber, the other one standard thermocouples type T use for ambient temperature measurement. The calibration was done in according to WI-T20 ( based on ASTM E145-94 ( Reapproved 2001 ) and AS2853-1986 ). All data show below were final values and the initial data from customer request. The temperature scale used was based on ITS - 90.

##### 2. Reference Standard Instrument :

Instrument	Model	Instrument No.	Certificate No.	Due Date
TC	TYPE T	TN161-TN170	T230773	10 April 2024
TC	TYPE T	TN171-TN180	T230773	10 April 2024
DATA LOGGER	34970A	T149	T230773	10 April 2024

##### 3. This certificate is traceable to :

National Institute of Metrology ( Thailand ) through Metrological Center ( NSC-TISI-TIS 17025 CALIBRATION 0244 )

##### 4. Condition of calibrated item : good

###### Equipment Description :

Time Constant 1 Hour 30 Minute At 3 °C

Fresh Air Damper ☐ Open ☐ Min ☐ Medium ☐ Max

☐ Close

☒ Not Available

##### 5. Adjustment :

( X ) without adjustment

( ) after adjustment

Approved By : 

FM-L15 118/18-08-66



## Metrology

SCI ECO Services Company Limited


33/2 Moo 3, T.Banpa, A.Kaengkhoh, Saraburi 18110, Thailand.



Certificate No. T232160

Page 3 of 4

### Calibration Report

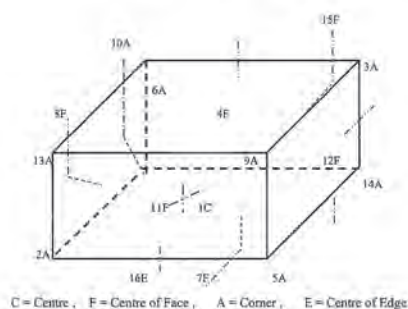

C = Centre , F = Centre of Face , A = Corner , E = Centre of Edge

1C = TN161	12F = TN172
2A = TN162	13A = TN173
3A = TN163	14A = TN174
4E = TN164	15F = TN175
5A = TN165	16E = TN176
6A = TN166	
7E = TN167	
8F = TN168	
9A = TN169	
10A = TN170	
11F = TN171	

Approved By : 

FM-L15 118/18-08-66





# Metrology SCI ECO Services Company Limited

33/2 Moo 3, T.Banpa, A.Kaengkhloi, Saraburi 18110, Thailand.



Certificate No. T232160

Page 4 of 4

## Calibration Report

### Measurement Results

Calibration Point	Average Standard Reading at each position (°C)									
	TN161	TN162	TN163	TN164	TN165	TN166	TN167	TN168	TN169	TN170
3.0	2.83	3.34	2.95	3.46	3.45	3.76	3.25	3.46	3.39	3.50
	TN171	TN172	TN173	TN174	TN175	TN176				
	3.33	3.39	3.15	3.43						

Chamber (Cooling Room)			Temperature Distribution				
Setting (°C)	Reading (°C)		Average (°C)	Stability (± °C)	Uniformity (°C)	Uncertainty (± °C)	Coverage Factor k
	Min	Max					
3.0	2.8	4.1	3.5	3.36	1.10	2.00	1.50

The calibration result apply only the above calibrated item.

The result of test was found accurate as shown on date and place of test only.

The reported expanded uncertainty is based on a standard uncertainty multiplied by a coverage factor k which for a t-distribution, providing a level of confidence of approximately 95 %.

Approved By:

*[Signature]*

TM-L1511818-08-66

REVIEW BY: *Auchawan S.*  
APPROVED BY: *Tanyatong M.*  
NEXT CAL DATE: 12 Jan 2025

### Certificate of Calibration

#### ICS-2100: Anion (ID#659)

This certificate is to verify that instrument below are calibrated

by Archemica Lab Co., Ltd.

ICS-2100 S/N: 15010977  
AS-HV S/N: 5450A36659

For

**ALS Laboratory Group (Thailand) Co., Ltd.**

Operator Signature: *Nutdanai* Date: Jan 12, 2024  
(Mr. Nutdanai Laekhwan)  
Application Chemist

Sartorius (Thailand) Co., Ltd.

120 Rama 9 Road, Huaykong, Huaykong, Bangkok 10250  
Tel: +66 2843 8301-6, e-mail: service.thailand@sartorius.com



REVIEW BY: *Jinda K.*  
APPROVED BY: *Shiruk P.*  
NEXT CAL DATE: 02/08/25

## Certificate of Calibration

Model Number: MSE224S-100-DU  
Description: Analytical Balance  
Serial Number: 0027405555  
ID No.: BXC\_EN0003  
Manufacturer: Sartorius

Certificate No.: 24BC0270  
Issued Date: Monday, August 05, 2024  
Reference No.: 240942  
Page No.: 1 of 2

Customer Name: ALS Laboratory Group (Thailand) Co., Ltd.  
104 Phatthanasarak 40/Phatthanasarak Rd, Khwaeng Suan Luang, Khet Suan Luang, Bangkok 10250

Calibrated Place: Lab Room

Calibrated By: Mr. Chuchit Jitthana

Calibration Date: Friday, August 02, 2024

Calibration Procedure No.: This calibration was conducted by  
Using in-house calibration procedure number (WI-003)  
Based on UKAS LAB 14 : 2019

### Metrological data:

Capacity: 220 g Repeatability: 0.0001 g

### Ambient Conditions:

Temperature: 23.0 °C ± 5.0 °C

Humidity: 55.0 % RH ± 10.0 % RH

Pressure: #

Equipment Condition: # Good Operate # Fail

Reasons for calibration: ☒ New installation ☒ Service / Repair ☒ Recalibration / Maintenance

Measurement Method: UKAS Publication Ref: Lab 14

The measurement uncertainty stated is the expanded uncertainty which is obtained from the standard uncertainty multiplied by the coverage factor (k=2) to provide a level of confidence of approximately 95%. It is determined in accordance with the Guide to Expression of Uncertainty in Measurement (GUM). The calibration certificate documents the traceability to National Standards, which realise the unit of measurement according to the International System of Units (SI). Report of Tolerance came from list of Sartorius Metrological Specifications.

### Traceability:

Model Number	Description	Traceability	Certificate No.	Due Date
YC5011-522-00	Sartorius weight set 100g - 5000g 12/YS011-522-00	YES	M23081975	23-Aug-2026
Testo 174 H	Thermo-Hygrometer, Testo 174H	ENTECH	11/17 661303/H861140	12-Nov-2024

This certificate entice and apply this equipment only.

This certificate may not be reproduced other than in full except with the prior written approval of the Verification Operation Division Sartorius (Thailand) Co., Ltd.

*[Signature]*  
Mr. Chuchit Jitthana (Technical Manager)



SOP FM 33 03 February 2022

Sartorius (Thailand) Co., Ltd.

120 Rama 9 Road, Huaykong, Huaykong, Bangkok 10250  
Tel: +66 2843 8301-6 Fax: +66 2843 8307, e-mail: service.thailand@sartorius.com

SARTORIUS

## Certificate of Calibration

Model Number: MSE224S-100-DU  
Description: Analytical Balance  
Serial Number: 0027405555  
ID No.: BXC\_EN0003  
Manufacturer: Sartorius

Certificate No.: 24BC0270  
Issued Date: Monday, August 05, 2024  
Reference No.: 240942  
Page No.: 2 of 2

### Calibration Results : Without Adjustment

Repeatability			Eccentricity (Off-center loading error)		
The repeatability is the ability of a weighing instrument to display nearly identical results under identical test conditions when the same load within a measurement range is placed repeatedly on the weighing pan in the same manner. The standard deviation is used to express repeatability uncertainty.			The off-center loading error is judged by the difference between the result of the load at 1/4 or 3/4 of maximum capacity placed in the middle of the weighing pan and between each of four additional measurement points (positions defined according to OIML R110).		
Nominal Value: (Low Load)	20.0000	200.0000	Nominal value:	100	g
20	20.0000	199.9999	Tolerance	0.0004	g
Tolerance:	20.0001	200.0000	Difference		
0.0001	20.0000	200.0000			
	20.0000	200.0000			
	20.0000	200.0000			
	20.0000	200.0000			
	20.0000	200.0000			
Nominal Value: (High Load)	200	200.0000			
200	200.0000	200.0000			
Tolerance:	20.0000	200.0000			
0.0001	20.0000	200.0000			
	20.0000	200.0000			
	20.0000	200.0000			
	20.0000	200.0000			
	20.0000	200.0000			
Standard Deviation	0.00004	0.00006			

Linearity				
The linearity, also called linearity error, describes the deviation of the characteristic curve of a weighing instrument from the linear slope.				
Tolerance	0.0002 g			
Nominal Value	Conventional Mass Value	Displayed Value	Deviation	Uncertainty
(g)	(g)	(g)	(g)	(g)
0.01	0.0100	0.0100	0.0000	0.00015
0.1	0.1000	0.1000	0.0000	0.00015
1	1.0000	1.0000	0.0000	0.00015
2	2.0000	2.0000	0.0000	0.00015
5	5.0000	5.0000	0.0000	0.00015
10	10.0000	10.0000	0.0000	0.00015
20	20.0000	20.0000	0.0000	0.00015
50	50.0000	50.0001	0.0001	0.00016
100	100.0000	100.0001	0.0001	0.00019
200	200.0000	200.0000	0.0000	0.00029
End of Report				

SOP FM 33 03 February 2022





## Metrology

SCI ECO Services Company Limited

33/2 Moo 3, T.Banpa, A.Kaengkhoh, Saraburi 18110, Thailand.  
Saraburi Tel : +66 3627 3096 Fax : +66 3627 3100  
Bangkok Tel : +668 9205 6851 , +669 8247 2360  
Website : www.scieco.co.th E-Mail : calibrate@scg.com



NIST-TIS-17025  
CALIBRATION 0244

Certificate No. T240904

Page 1 of 3

### Certificate of Calibration

Equipment : Chamber ( Oven )

Manufacturer : Memmert

Model : UF 450

Serial No. : B717.0531

Customer Code : BKK\_EN0273

ID No. : T8042A4

Customer : ALS Laboratory Group (Thailand) Co.,Ltd.

104 Phatthanakan 40, Phatthanakan Rd., Khwaeng Phatthanakan,

Khet Suan Luang, Bangkok 10250


Customer Location : Laboratory (Oven Room)

Date of Receipt : 08 May 2024

Calibrated By : Preecha Phisassuthikul ( Temperature Calibration Manager )

Approved By :                      / Nuafun Sungsom (Metrology Manager)

Date of Issue : 23 MAY 2024

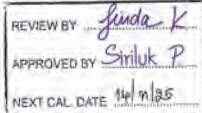

The uncertainties are for a confidence probability of approximately 95%.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standard laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the Metrology.

FM-L14 11/18-08-66



## Metrology

SCI ECO Services Company Limited

33/2 Moo 3, T.Banpa, A.Kaengkhoh, Saraburi 18110, Thailand.



NIST-TIS-17025  
CALIBRATION 0244

Certificate No. T240904

Page 2 of 3

### Calibration Report

Equipment : Chamber ( Oven )

Date of Calibration : 14 May 2024

Environment : Temperature : 26.5-28.1 °C

Line Voltage : 226.7-229.8 V

Relative Humidity : 51 - 57 %RH

#### Condition of this results of calibration :

1. This equipment was calibrated by insert nine resistance thermometer detectors into its chamber , the other one resistance thermometer detector use for ambient temperature measurement . The calibration was done in according to WI-T20 ( based on ASTM E145-94 ( Reapproved 2001) and AS2853-1986 ) . All data show below were final values and the initial data from customer request . The temperature scale used was based on ITS - 90 .

#### 2. Reference Standard Instrument :

Instrument	Model	Instrument No.	Certificate No.	Due Date
RTD	100 ohm	21-(CH1-10)	T231955	17 November 2024
DATA LOGGER	34970A	T121	T231955	17 November 2024

#### 3. This certificate is traceable to :

National Institute of Metrology ( Thailand ) through Metrological Center ( NSC-TIS-TIS 17025 CALIBRATION 0244 ) .

#### 4. Condition of calibrated item : good

##### Equipment Description :

Time Constant	1 Hour	30 Minute	At 104 °C
Fresh Air Damper	<input type="checkbox"/> Open <input type="checkbox"/> Min <input type="checkbox"/> Medium <input type="checkbox"/> Max		
	<input type="checkbox"/> Close		
	<input checked="" type="checkbox"/> Not Available		

#### 5. Adjustment :

( X ) without adjustment ( ) after adjustment

Approved By :                     

FM-L15 11/18-08-66



## Metrology

SCI ECO Services Company Limited

33/2 Moo 3, T.Banpa, A.Kaengkhoh, Saraburi 18110, Thailand.




NIST-TIS-17025  
CALIBRATION 0244

Certificate No. T240904

Page 3 of 3

### Calibration Report

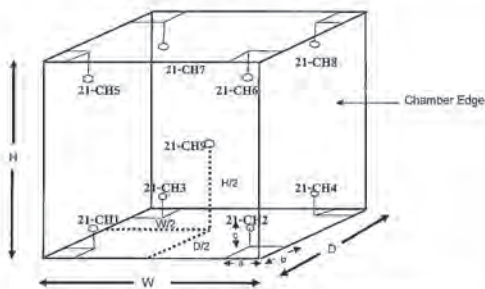

#### Remark :

Internal Dimensions of Chamber : W (Width) = 104 cm , H (Height) = 72 cm and D (Depth) = 60 cm .  
Size of Installed Standard sensor number 21-CH1 to number 21-CH11 : a = 5 cm , b = 5 cm and c = 5 cm .  
Size of Installed Standard sensor number 21-CH9 : W/2 = 104 cm/2 , H/2 = 72 cm/2 and D/2 = 60 cm/2

#### Measurement Results

Average Standard Reading at each position ( °C )								
Calibration Point	21-CH1	21-CH2	21-CH3	21-CH4	21-CH5	21-CH6	21-CH7	21-CH8
104	103.4	105.0	103.7	103.6	103.3	104.6	103.3	104.0
180	179.5	181.1	179.2	179.5	179.0	181.3	179.3	179.9

Chamber ( Oven )		Temperature Distribution					
Setting ( °C )	Reading ( °C )		Average ( °C )	Stability ( °C )	Uniformity ( °C )	Uncertainty ( °C )	Coverage Factor k
	Min	Max					
104.0	103.9	104	104.0	0.14	1.27	0.84	2.00
180.0	179.9	180.1	180.0	0.39	2.29	0.76	2.00

\* The quoted uncertainty excludes "uniformity"

The calibration result apply only the above calibrated item.

The result of test was found accurate as shown on date and place of test only.

The reported expanded uncertainty is based on a standard uncertainty multiplied by a coverage factor k which for a t-distribution, providing a level of confidence of approximately 95 % .

End of Certificate

Approved By :                     

FM-L15 11/18-08-66



## Metrological Center

SCI ECO Services Company Limited

33/2 Moo 3, T.Banpa, A.Kaengkhoh, Saraburi 18110

Telephone : +66 2 586 5792-4 Fax : +66 2 586 5109

Website : www.scieco.co.th E-Mail : calibrate@scg.co.th

Certificate No. T240742

Page 1 of 5

### Certificate of Calibration

Equipment : Digestion Unit

Manufacturer : SCP Science

Model : DigiPRER HT

Serial No. : HTC1120480658

Customer Code : BKK\_EN0366

ID No. : T2635A5

Customer : ALS Laboratory Group (Thailand) Co.,Ltd.

104 Phatthanakan 40, Phatthanakan Rd., Khwaeng Phatthanakan,

Khet Suan Luang, Bangkok 10250


Customer Location : Wet Chemistry Lab 1

Date of Receipt : 11 April 2024

Calibrated By : Sujjar Naknakred ( Site Calibration Manager )

Approved By :                      / Boonchai Suriyawong ( Site Calibration Manager )

Date of Issue : 02 MAY 2024

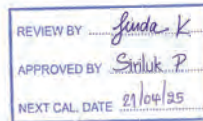

The uncertainties are for a confidence probability of approximately 95%.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standard laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the Metrological Center.

FM-L12 10/30-05-57





## Metrological Center

### SCI ECO Services Company Limited

33/2 Moo 3, T.Banpa, A.Kaengkhoh, Saraburi 18110  
Telephone : +66 2 586 5792-4 Fax : +66 2 586 5109  
Website : www.scieco.co.th E-Mail : calibrate@scg.co.th

Certificate No. T240742

Page 2 of 5

## Calibration Report

Equipment : Digestion Unit  
Date of Calibration : 21 April 2024  
Environment : Temperature : 23.9 - 26.3 °C  
Line Voltage : 221.8 - 225.9 V  
Relative Humidity : 55 - 65 %RH

### Condition of this results of calibration :

- This equipment was calibrated by insert four standard thermocouples type S into its chamber, the other one thermocouple type T use for ambient temperature measurement, The calibration was done in according to WI-T10.  
was based on ITS - 90.
- Reference Standard Instrument :

Instrument	Model	Instrument No.	Certificate No.	Due Date
TC	Type S	M20A2-(CH11-CH14)	T230886	09 May 2024
DATA LOGGER	34970A	T47	T230886	09 May 2024
- This certificate is traceable to :  
National Institute of Metrology ( Thailand ) through Metrological Center ( NSC-TISI-TIS 17025 CALIBRATION 0244 )
- Condition of calibrated item : good  
Equipment Description :  
Time Constant ☒ 1 Hour ☒ 6 Minute At 380 °C  
Fresh Air Damper ☐ Open ☐ Min ☐ Medium ☐ Max  
☐ Close  
☒ Not Available
- Adjustment :  
( X ) without adjustment ( ) after adjustment

Approved By:

FM-L13 108/30-05-57



## Metrological Center


### SCI ECO Services Company Limited

33/2 Moo 3, T.Banpa, A.Kaengkhoh, Saraburi 18110  
Telephone : +66 2 586 5792-4 Fax : +66 2 586 5109  
Website : www.scieco.co.th E-Mail : calibrate@scg.co.th

Certificate No. T240742

Page 3 of 5

## Calibration Report

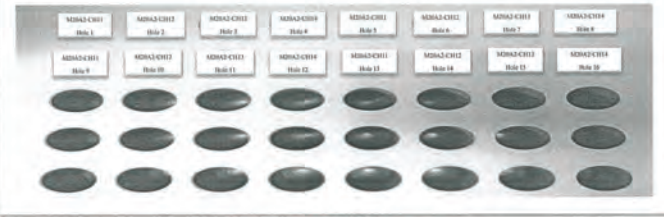

FRONT

### Measurement Results

Cal. Point	Setting	Reading	STD.	Position of Standards at Block							
(°C)	(°C)	(°C)	Reading	M20A2-CH11 Block 1	M20A2-CH12 Block 2	M20A2-CH13 Block 3	M20A2-CH14 Block 4	M20A2-CH15 Block 5	M20A2-CH16 Block 6	M20A2-CH17 Block 7	M20A2-CH18 Block 8
380.0	380.0	379.2 - 380.5	Max °C	378.7	378.9	377.9	378.7	380.5	379.8	378.7	377.4
			Min °C	378.2	378.5	377.5	378.2	380.1	379.3	378.3	376.9
			Average °C	378.4	378.7	377.7	378.4	380.3	379.6	378.5	377.2
			Stability ± °C	0.2	0.2	0.2	0.2	0.2	0.3	0.2	0.2

Cal. Point	Setting	Reading	STD.	Position of Standards at Block							
(°C)	(°C)	(°C)	Reading	M20A2-CH11 Block 9	M20A2-CH12 Block 10	M20A2-CH13 Block 11	M20A2-CH14 Block 12	M20A2-CH15 Block 13	M20A2-CH16 Block 14	M20A2-CH17 Block 15	M20A2-CH18 Block 16
380.0	380.0	379.2 - 380.5	Max °C	378.4	378.6	379.2	379.6	381.9	380.6	379.1	378.1
			Min °C	377.8	378.2	378.7	379.2	381.4	379.9	378.3	377.2
			Average °C	378.1	378.4	379.6	379.4	381.6	380.3	378.7	377.7
			Stability ± °C	0.3	0.2	0.2	0.2	0.3	0.4	0.4	0.5

Approved By:

FM-L13 108/30-05-57



## Metrological Center

### SCI ECO Services Company Limited

33/2 Moo 3, T.Banpa, A.Kaengkhoh, Saraburi 18110  
Telephone : +66 2 586 5792-4 Fax : +66 2 586 5109  
Website : www.scieco.co.th E-Mail : calibrate@scg.co.th

Certificate No. T240742

Page 4 of 5

## Calibration Report



FRONT

### Measurement Results

Cal. Point	Setting	Reading	STD.	Position of Standards at Block							
(°C)	(°C)	(°C)	Reading	M20A2-CH11 Block 17	M20A2-CH12 Block 18	M20A2-CH13 Block 19	M20A2-CH14 Block 20	M20A2-CH15 Block 21	M20A2-CH16 Block 22	M20A2-CH17 Block 23	M20A2-CH18 Block 24
380.0	380.0	379.2 - 380.5	Max °C	378.9	379.2	379.5	380.1	382.1	381.0	378.9	377.8
			Min °C	378.2	378.6	379.1	379.6	381.7	380.2	378.3	377.2
			Average °C	378.5	378.9	379.3	379.8	381.9	380.6	378.6	377.5
			Stability ± °C	0.3	0.3	0.2	0.2	0.2	0.4	0.3	0.3

Cal. Point	Setting	Reading	STD.	Position of Standards at Block							
(°C)	(°C)	(°C)	Reading	M20A2-CH11 Block 25	M20A2-CH12 Block 26	M20A2-CH13 Block 27	M20A2-CH14 Block 28	M20A2-CH15 Block 29	M20A2-CH16 Block 30	M20A2-CH17 Block 31	M20A2-CH18 Block 32
380.0	380.0	379.2 - 380.5	Max °C	378.5	378.1	378.0	378.6	380.7	379.7	377.7	380.9
			Min °C	378.2	377.8	377.7	378.1	380.3	379.0	377.2	380.4
			Average °C	378.4	378.0	377.8	378.4	380.5	379.4	377.5	380.6
			Stability ± °C	0.2	0.2	0.2	0.2	0.2	0.2	0.2	0.3

Approved By:

FM-L13 108/30-05-57



## Metrological Center


### SCI ECO Services Company Limited

33/2 Moo 3, T.Banpa, A.Kaengkhoh, Saraburi 18110  
Telephone : +66 2 586 5792-4 Fax : +66 2 586 5109  
Website : www.scieco.co.th E-Mail : calibrate@scg.co.th

Certificate No. T240742

Page 5 of 5

## Calibration Report

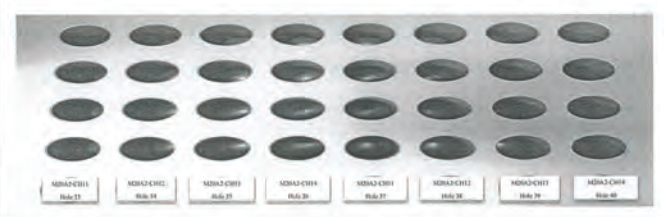

FRONT

### Measurement Results

Cal. Point	Setting	Reading	STD.	Position of Standards at Block							
(°C)	(°C)	(°C)	Reading	M20A2-CH11 Block 33	M20A2-CH12 Block 34	M20A2-CH13 Block 35	M20A2-CH14 Block 36	M20A2-CH15 Block 37	M20A2-CH16 Block 38	M20A2-CH17 Block 39	M20A2-CH18 Block 40
380.0	380.0	379.2 - 380.5	Max °C	378.3	377.9	378.7	379.5	381.6	380.5	378.4	378.0
			Min °C	378.0	377.6	378.4	379.1	381.2	380.0	378.1	377.6
			Average °C	378.2	377.8	378.6	379.3	381.4	380.3	378.2	377.8
			Stability ± °C	0.2	0.2	0.2	0.2	0.2	0.2	0.2	0.2

The expanded uncertainty of temperature measurement was  $\pm 1.87$  °C

The calibration result apply only the above calibrated item.

The result of test was found accurate as shown on date and place of test only.

The reported expanded uncertainty is based on a standard uncertainty multiplied by a coverage factor k=2, providing a level of confidence of approximately 95 %.

Approved By:

FM-L13 108/30-05-57





Calibration Results:

Before Adjustment				
Standard	Unit Under Calibration	Reading	Correction	Coverage Factor
Conductivity Solution	Reading	Correction	Coverage Factor	Uncertainty (±)
84,000 µS/cm	92.64 µS/cm	-8.640 µS/cm	2.00	0.68 µS/cm
1413.0 µS/cm	1423 µS/cm	+10.0 µS/cm	2.00	11 µS/cm
12,880 mS/cm	12.81 mS/cm	0.070 mS/cm	2.00	0.10 mS/cm

After Adjustment : at 84.0 µS/cm, 1413 µS/cm, 12.88 mS/cm

Standard	Unit Under Calibration	Reading	Correction	Coverage Factor
Conductivity Solution	Reading	Correction	Coverage Factor	Uncertainty (±)
84,000 µS/cm	84.03 µS/cm	-0.030 µS/cm	2.00	0.68 µS/cm
1413.0 µS/cm	1414 µS/cm	-1.0 µS/cm	2.00	11 µS/cm
12,880 mS/cm	12.88 mS/cm	0.020 mS/cm	2.00	0.098 mS/cm

The End of Certificate

DKSH Technology Limited  
2533 Sukhumvit Road, Bangkok, Thailand 10260  
Phone: +66 2639 7100 Email: info@dksh.com

Delivering Growth - In Asia and Beyond.

CAL-EM-023-012/17.640.2022

ใบตรวจสอบสภาพเครื่องวัดสิ่งแวดล้อม

เลขที่ใบงาน: WO-00012682

ชนิดเครื่องมือ: CONDUCTIVITY METER

รุ่น: ORION STAR A215

หมายเลขเครื่อง: X58031

ตรวจสอบ (วัน)		รายการตรวจเช็ค	ตรวจสอบ (ผู้)		หมายเหตุ
25 Dec 2023			25 Dec 2023		
ปกติ	ไม่ปกติ		ปกติ	ไม่ปกติ	
General					
<input checked="" type="checkbox"/>	<input type="checkbox"/>	1. ความสมบูรณ์เครื่อง	<input checked="" type="checkbox"/>	<input type="checkbox"/>	
<input checked="" type="checkbox"/>	<input type="checkbox"/>	2. ความสะอาด (ช่องใส่ตัวอย่าง, ภายใน-นอกเครื่อง)	<input checked="" type="checkbox"/>	<input type="checkbox"/>	
<input checked="" type="checkbox"/>	<input type="checkbox"/>	3. สวิตช์ ปิด-เปิด เครื่อง (On-Off Switch)	<input checked="" type="checkbox"/>	<input type="checkbox"/>	
<input checked="" type="checkbox"/>	<input type="checkbox"/>	4. ปุ่มกด (Keypad)	<input checked="" type="checkbox"/>	<input type="checkbox"/>	
<input checked="" type="checkbox"/>	<input type="checkbox"/>	5. หน้าจอ (Display, Screen Contrast)	<input checked="" type="checkbox"/>	<input type="checkbox"/>	
Spectrophotometer					
<input type="checkbox"/>	<input type="checkbox"/>	6. แบตเตอรี่สำรอง (Battery Backup) >= 2.5 VDC	<input type="checkbox"/>	<input type="checkbox"/>	
<input type="checkbox"/>	<input type="checkbox"/>	7. ควบคุมแสงความยาวคลื่น (Wavelength Control)	<input type="checkbox"/>	<input type="checkbox"/>	
<input type="checkbox"/>	<input type="checkbox"/>	8. ความยาวคลื่น (Wavelength Check)	<input type="checkbox"/>	<input type="checkbox"/>	
<input type="checkbox"/>	<input type="checkbox"/>	9. แหล่งกำเนิดแสง (UV < 3,000 hour)	<input type="checkbox"/>	<input type="checkbox"/>	
<input type="checkbox"/>	<input type="checkbox"/>	10. แหล่งกำเนิดแสง (Visible < 5,000 hour)	<input type="checkbox"/>	<input type="checkbox"/>	
<input type="checkbox"/>	<input type="checkbox"/>	11. ช่องใส่หลอดตัวอย่าง (Carousel Module)	<input type="checkbox"/>	<input type="checkbox"/>	
pH Meter and Conductivity Meter					
<input type="checkbox"/>	<input type="checkbox"/>	12. บัสบาร์ (Electrode and Connection Cable)	<input type="checkbox"/>	<input type="checkbox"/>	
<input type="checkbox"/>	<input type="checkbox"/>	13. ระดับสารละลายใน Electrode (Level KCl)	<input type="checkbox"/>	<input type="checkbox"/>	
<input type="checkbox"/>	<input type="checkbox"/>	14. ฝาปิดกันน้ำ Electrode (Dust Protection Hood)	<input type="checkbox"/>	<input type="checkbox"/>	
<input type="checkbox"/>	<input type="checkbox"/>	15. ขาตั้งยึดเครื่อง (Stand)	<input type="checkbox"/>	<input type="checkbox"/>	
Turbidimeter					
<input type="checkbox"/>	<input type="checkbox"/>	16. ค่าความขุ่นที่ต่ำสุด (No Sample)	<input type="checkbox"/>	<input type="checkbox"/>	
<input type="checkbox"/>	<input type="checkbox"/>	17. ระดับการสั่นสะเทือนเครื่อง (>= 2.5 ไม่นับ 3.0)	<input type="checkbox"/>	<input type="checkbox"/>	
Automatic titrator					
<input type="checkbox"/>	<input type="checkbox"/>	18. สลัก Piston Burettes	<input type="checkbox"/>	<input type="checkbox"/>	
<input type="checkbox"/>	<input type="checkbox"/>	19. Function Rinsing and Dosing	<input type="checkbox"/>	<input type="checkbox"/>	
<input type="checkbox"/>	<input type="checkbox"/>	20. ระบบท่อจ่ายและดูดจากถังพัก	<input type="checkbox"/>	<input type="checkbox"/>	

ผู้ตรวจ :

Mr.Sivaporn Srijan  
Service Engineer

DKSH Technology Limited  
2533 Sukhumvit Road, Bangkok, Thailand 10260  
Phone: +66 2639 7100 Email: info@dksh.com

Delivering Growth - In Asia and Beyond.

CAL-FABR14-03/20 Jul.2022




TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES  
534/4 PATTANAKARN ROAD SOI 16, SUANLUANG, SUANLUANG BANGKOK 10250  
TEL: 0-2717-3000 FAX: 0-2719-9488

Cert.No.: 23TW243  
Page: 1 of 2

Certificate of Testing

Equipment : DO Meter  
Manufacturer : YSI  
Model : 5000-230V  
Serial No. : 08J101147  
ID No. : BKK\_EN0017  
Received Date : 15 November 2023  
Test Date : 16 November 2023  
Reference : 2311-0505DSC-4  
Submitted by : ALS Laboratory Group (Thailand) Co., Ltd.  
104 Phatthanakan 40, Phatthanakan Rd.,  
Khuwaeng Phatthanakan, Khet Suan Luang,  
Bangkok 10250 Thailand  
Laboratory Condition : Temperature (25 ± 5) °C  
Humidity (50 ± 20) %  
Test Procedure : In-house method : CP-CH9  
by Comparison Technique with Azide Modification Method  
Tested by : Waijaiak Srinhean  
Approved by :   
Approved Signatory  
(✓) Saitrip Meangmai  
( ) Warakorn Lamgagatrakul  
( ) Ponpan Paipim  
Issue Date : 17 November 2023

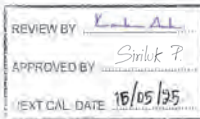

Cert.No.: 23TW243  
Page: 2 of 2

Condition of this result of calibration

1. Reference Standard Instruments :  
This certification is traceable to the International System of Unit through the reference standards  
laboratory of Industrial Calibration Center, Technology Promotion Association (Thailand-Japan).

Instruments	Serial No.	ID No.	Certificate No.	Due Date
1) Burette	-	1308J10	23CG1172	22 Mar 2025
2) Balance	1124013362	140RC006	23MM16	20 Feb 2024
Standard Material :-				
Material	Manufacturer	Lot No.	Assay	
Sodium Thiosulfate pentahydrate	Merck	AM1763316	100.2%	

Result : Dissolved Oxygen Meter Adjustment With Air 100 %  
Dissolved Oxygen Probe No.: 16K100496

Titration Method (Azide Modification Method) (mg/L)	DO Meter Reading (mg/L)	Standard Deviation (mg/L)
8.18	8.18	0.0055

This report was certified only for the instrument we tested. It is allowable to use for study  
the system efficiency. The environmental impact control and present to organization it may concerned  
Intend to use for advertising and referral purpose is prohibited. This report may not be reproduced  
other in full without written approval of the laboratory

-080-

1190297

0328589





Cert. No.: 23LM192  
Page.: 1 of 2

## Certificate of Calibration

Equipment : DO Meter with Sensor  
Manufacturer : YSI  
Model : 5000-230V  
Serial No. : 09J101147  
ID No. : BKK\_EN0017  
Submitted by : ALS Laboratory Group (Thailand) Co., Ltd.  
104 Phatthanakan 40, Phatthanakan Rd.,  
Khwaeng Phatthanakan, Khet Suan Luang,  
Bangkok 10250 Thailand  
Location : TPA Chemistry Calibration Laboratory  
Received Order : 15 November 2023  
Calibrated Date : 16 November 2023  
Ambient Temperature :  $(28 \pm 10) ^\circ\text{C}$   
Relative Humidity :  $(50 \pm 30) \%$   
AC Line Voltage :  $(220 \pm 22) \text{ V}$   
Calibrated by : Kunchit Prompru  
Approved by :   
( ) Pomthippa Tamayakul  
( ) Ponpan Palpin  
(✓) Suwit Imjai  
Issue Date : 17 November 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full except with the prior written approval of the Association of Corporate Services (ACS) - Equipment Calibration and Testing Services.

A 0060730



Equipment : DO Meter with Sensor  
Condition As-Received : Used Item  
Reference : 2311-0505DSC-10  
Procedure Used :-

Cert. No.: 23LM192  
Page.: 2 of 2

Calibration were conducted using in-house calibration procedure CP-OT01 according to comparison with Industrial Platinum Resistance Thermometer (IPRT) into Temperature Bath.  
The temperature scale used was based on ITS-90.

### Condition of this result of calibration

#### 1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Digital Thermometer	3240076	231305	TPA	15 Mar 2024

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (°) Without Adjustment

Function : Temperature measurement.

This instrument was connected with temperature sensor, S/N : 16K100496

Calibration Point (°C)	Immersion Depth (mm)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty (± °C)	Coverage Factor k
20.0	60	19.997	19.93	-0.067	0.15	2.00

UUC\* : Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor  $k$ , providing a level of confidence of approximately 95%.

-000-

a 1190298



## QUALITY CALIBRATION CO.,LTD.

235 Petchkasem 63/2 Road, Laksoeng, Bangkok, Bangkok 10160  
Tel (662) 421-5402, (662) 444-0152-3, Fax (662) 809-4584  
www.qcalibration.com



CERTIFICATE No : 24T2852  
REFERENCE No : 72619-8

PAGE : 1 OF 2

## Certificate of Calibration

EQUIPMENT : COOLED INCUBATOR  
MANUFACTURER : MEMMERT  
MODEL : ICP750  
SERIAL No : F819.0021  
ID No : BKK\_EN0304  
CONDITION AS RECEIVED : USED ITEM  
SUBMITTED BY : ALS LABORATORY GROUP (THAILAND) CO., LTD.  
104 PHATTHANAKAN 40, PHATTHANAKAN  
RD., KHWAENG PHATTHANAKAN, KHET SUAN  
LUANG, BANGKOK 10250, THAILAND  
CALIBRATED BY : CHAICHARN CH.  
CALIBRATION DATE : 20-Mar-24  
APPROVED BY :   
PONGSAK J.  
ISSUED DATE : 21-Mar-24  
RECEIVED DATE : 20-Mar-24

REVIEW BY :   
APPROVED BY :   
NEXT CAL DATE : 30/03/25



## QUALITY CALIBRATION CO.,LTD.

235 Petchkasem 63/2 Road, Laksoeng, Bangkok, Bangkok 10160  
Tel (662) 421-5402, (662) 444-0152-3, Fax (662) 809-4584

CERTIFICATE No : 24T2852

PAGE : 2 OF 2

## Calibration Report

EQUIPMENT : COOLED INCUBATOR  
MANUFACTURER : MEMMERT  
MODEL : ICP750  
ID No : BKK\_EN0304  
RECEIVED DATE : 20-Mar-24  
AMBIENT TEMPERATURE :  $26 ^\circ\text{C} \pm 1 ^\circ\text{C}$   
S/N : F819.0021  
CALIBRATION DATE : 20-Mar-24  
RELATIVE HUMIDITY :  $54 \text{ \%RH} \pm 10 \text{ \%RH}$

### CONDITION OF THIS RESULTS OF CALIBRATION


1. THIS INSTRUMENT WAS CALIBRATED ACCORDING TO TLA5 G-20 BY COMPARISON WITH CALIBRATED THERMOCOUPLE TYPE K UNDER NO LOAD CONDITION. THE THERMOCOUPLES WERE PLACED ON NINE POINTS AND LOCATED ONE THERMOCOUPLE IN EACH OF THE EIGHT CORNERS OF THE CHAMBER AND WAS AWAY FROM THE EACH WALL OF 5 cm TO 10 cm. AND PLACED THE NINTH THERMOCOUPLE WITHIN 2.5 cm. OF THE GEOMETRIC CENTER OF THE CHAMBER. THE UNIFORMITY WAS MEASURED BETWEEN REFERENCE PROBE AND OTHER PROBES AT THE SAME TIME.

#### 2. REFERENCE STANDARD INSTRUMENTS :-

INSTRUMENT	MODEL	SERIAL No	CERTIFICATE No	DUE DATE
1) DATA LOGGER WITH TC TYPE K	HYDRA 2635A	7286308	2316641	14-Jul-24

3. THE CERTIFICATE IS VALID FOR THE ITEM CALIBRATED AS SHOWN ON THE DATE AND PLACE OF CALIBRATION ONLY.  
4. THIS RESULT EXCLUDE LONG TERM STABILITY OF THE UNIT UNDER CALIBRATION.  
5. THIS CERTIFICATE IS TRACEABLE TO THE INTERNATIONAL SYSTEM OF UNIT MAINTAINED AT:-  
- NATIONAL INSTITUTE OF METROLOGY (THAILAND) THROUGH QUALITY CALIBRATION CO., LTD.

### RESULT OF CALIBRATION :- WITHOUT ADJUSTMENT

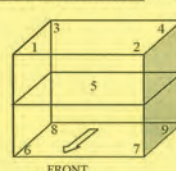

GENERAL INFORMATION  
Overall Ambient Temperature around the Chamber (°C) variation : 1  
Overall Line Voltage (V) variation : 5  
Instrument Condition : Normal

#### CHAMBER PERFORMANCE

Controller Temperature (°C)	Indicating Temperature (°C)	Temperature Stability (±°C)	Temperature Uniformity (°C)	Overall Variation (°C)
20.0	20.0	0.16	0.21	0.41

#### TEMPERATURE MEASUREMENT ACCURACY TEST

Controller	Indicating	Measured Temperature (°C) at Spread Locations									Uncertainty (±°C)
Temp (°C)	Temp (°C)	#1	#2	#3	#4	Ref. 5	#6	#7	#8	#9	
20.0	20.0	19.88	19.93	19.87	19.86	19.98	19.94	19.94	19.89	19.91	0.42

NOTE 1 : THE UNCERTAINTY OF MEASUREMENT EXCLUDED TEMPERATURE UNIFORMITY OF THE CHAMBER.  
NOTE 2 : LOCATION 5 WAS REFERENCE LOCATION.

NOTE 3 : THIS CALIBRATION WAS CARRIED OUT AT THE CUSTOMER'S PLACE AT LABORATORY AREA.  
THE REPORTED UNCERTAINTY OF MEASUREMENT WAS BASED ON A STANDARD UNCERTAINTY MULTIPLIED BY A COVERAGE FACTOR  $k=2$ , PROVIDING A LEVEL OF CONFIDENCE APPROXIMATELY 95%.

END OF CALIBRATION REPORT






Certificate No. T240741

Page 1 of 5

### Certificate of Calibration

Equipment : HOT BLOCK  
Manufacturer : Environmental Express  
Model : B3000-240  
Serial No. : 2017CODW116  
Customer Code : BKK\_EN0222  
ID No. : T6769A4  
Customer : ALS Laboratory Group (Thailand) Co.,Ltd.  
104 Phatthanakan 40, Phatthanakan Rd., Khwaeng Phatthanakan,  
Khet Suan Luang, Bangkok 10250  
Customer Location : Wet Chemistry Lab2  
Date of Receipt : 11 April 2024  
Calibrated By : Sane Musikanwan (Site Calibration Manager)  
Approved By :  / Sujar Naknakred (Site Calibration Manager)  
Date of Issue : 23 APR 2024

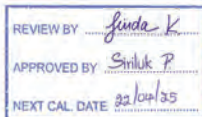

The uncertainties are for a confidence probability of approximately 95%.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standard laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the Metrological Center.

FM-L12 108/30-05-57



Certificate No. T240741

Page 2 of 5

### Calibration Report

Equipment : HOT BLOCK  
Date of Calibration : 22 April 2024  
Environment : Temperature : 22.9-24.4 °C  
Line Voltage : 222.7-227.8 V  
Relative Humidity : 55 - 65 %RH

#### Condition of this results of calibration :

- This equipment was calibrated by insert 20 standard thermocouples type T into its chamber, the other one standard thermocouples type T use for ambient temperature measurement. The calibration was done in according to WI-T20 (based on ASTM E145-94 (Reapproved 2001) and AS2853-1986). All data show below were final values and the initial data from customer request. The temperature scale used was based on ITS - 90.
- Reference Standard Instrument :

Instrument	Model	Instrument No.	Certificate No.	Due Date
TC	TYPE T	TN21-TN30	T240235	17 February 2025
TC	TYPE T	TN31-TN40	T240235	17 February 2025
DATA LOGGER	34970A	T195	T240235	17 February 2025
- This certificate is traceable to :  
National Institute of Metrology ( Thailand ) through Metrological Center ( NSC-TISI-TIS 17025 CALIBRATION 0244 )
- Condition of calibrated item : good  
Equipment Description :  
Time Constant 1 Hour 10 Minute At 150 °C  
Fresh Air Damper ☐ Open ☐ Min ☐ Medium ☐ Max  
☐ Close  
☒ Not Available
- Adjustment :  
( X ) without adjustment ( ) after adjustment

Approved By. 


FM-L13 108/30-05-57



Certificate No. T240741

Page 3 of 5

### Calibration Report

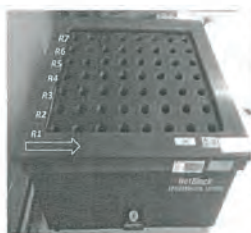

Row	Hole							
R7	H49	H50	H51	H52	H53	H54	H55	H56
R6	H41	H42	H43	H44	H45	H46	H47	H48
R5	H33	H34	H35	H36	H37	H38	H39	H40
R4	H25	H26	H27	H28	H29	H30	H31	H32
R3	H17	H18	H19	H20	H21	H22	H23	H24
R2	H9	H10	H11	H12	H13	H14	H15	H16
R1	H1	H2	H3	H4	H5	H6	H7	H8

II: STANDARD THERMOCOUPLE TYPE T

H1	=	TN21	H9	=	TN29	H17	=	TN37	H25	=	TN25	H33	=	TN33	H41	=	TN21	H49	=	TN29
H2	=	TN22	H10	=	TN30	H18	=	TN38	H26	=	TN26	H34	=	TN34	H42	=	TN22	H50	=	TN30
H3	=	TN23	H11	=	TN31	H19	=	TN39	H27	=	TN27	H35	=	TN35	H43	=	TN23	H51	=	TN31
H4	=	TN24	H12	=	TN32	H20	=	TN40	H28	=	TN28	H36	=	TN36	H44	=	TN24	H52	=	TN32
H5	=	TN25	H13	=	TN33	H21	=	TN21	H29	=	TN29	H37	=	TN37	H45	=	TN25	H53	=	TN33
H6	=	TN26	H14	=	TN34	H22	=	TN22	H30	=	TN30	H38	=	TN38	H46	=	TN26	H54	=	TN34
H7	=	TN27	H15	=	TN35	H23	=	TN23	H31	=	TN31	H39	=	TN39	H47	=	TN27	H55	=	TN35
H8	=	TN28	H16	=	TN36	H24	=	TN24	H32	=	TN32	H40	=	TN40	H48	=	TN28	H56	=	TN36

Approved By. 

FM-L13 108/30-05-57



Certificate No. T240741

Page 4 of 5

### Calibration Report

#### Measurement Results

		Average Standard Reading at each position (°C)									
Calibration Point		TN21	TN22	TN23	TN24	TN25	TN26	TN27	TN28	TN29	TN30
150	Setting	Max	150.38	149.73	150.12	150.36	150.56	150.00	150.48	150.25	150.56
	Min	150.43	149.47	149.87	150.16	150.31	149.77	150.25	150.02	149.41	149.41
	Average	150.23	149.59	149.96	150.24	150.41	149.87	150.36	150.12	150.45	149.51
	Setting	Max	150.17	150.28	150.28	150.37	150.09	149.96	149.86	149.75	150.63
	Min	149.94	150.03	150.01	150.18	149.88	149.69	149.68	149.57	150.41	149.96
	Average	150.04	150.14	150.13	150.27	149.98	149.81	149.77	149.63	150.51	150.03
	Setting	Max	150.28	150.18	149.87	149.57	150.18	149.90	150.59	149.66	150.39
	Min	150.00	149.94	149.67	149.39	149.88	149.38	150.32	149.34	150.11	149.84
	Average	150.14	150.07	149.77	149.49	150.04	149.75	150.48	149.52	150.26	149.97
	Setting	Max	150.38	149.71	150.18	149.97	150.03	150.05	150.21	150.07	150.02
	Min	150.12	149.49	149.87	149.66	149.71	149.71	149.89	149.79	149.76	149.73
	Average	150.26	149.61	150.04	149.82	149.90	149.89	150.05	149.94	149.91	149.84
	Setting	Max	150.37	150.20	150.20	150.44	150.67	149.85	150.31	149.90	150.36
	Min	150.11	149.99	150.04	150.26	150.49	149.69	150.12	149.78	150.20	149.40
	Average	150.25	150.12	150.14	150.34	150.57	149.78	150.20	149.83	150.29	149.52
	Setting	Max	150.18	150.02	149.95	150.26	149.92	149.69			
	Min	150.06	149.88	149.79	150.12	149.80	149.58				
	Average	150.13	149.95	149.89	150.18	149.84	149.64				

Approved By. 

FM-L13 108/30-05-57





## Metrological Center

SCI ECO Services Company Limited

33/2 Moo 3, T.Banpa, A.Kaengkhoi, Saraburi 18110  
Telephone : +66 2 586 5792-4 Fax : +66 2 586 5109  
Website : www.scieco.co.th E-Mail : calibrate@scg.co.th

Certificate No. T240741

Page 5 of 5

### Calibration Report

#### Measurement Results

HOT BLOCK		Temperature Distribution		
Setting (°C)	Reading (°C)		Stability (± °C)	Uncertainty (± °C)
	Min, Max	Average		
150.0	150, 150.1	150.0	0.20	0.83

The calibration result apply only the above calibrated item.

The result of test was found accurate as shown on date and place of test only.

The reported expanded uncertainty is based on a standard uncertainty multiplied by a coverage factor k=2, providing a level of confidence of approximately 95 %.

Approved By

FM-L13 108-30-05-57



Bara Scientific Co., Ltd.  
968 U Chu Liang Building Floor7 Ramat Road  
Silom Bangkok Bangkok Thailand 10500  
Tel : 02-6324300 Fax : 02-6375496-7  
www.barascientific.com



## Certificate of Calibration

Number of Page(s) 1 of 3

Certificate No. BSCC-UV-374/24  
Equipment UV/Vis Spectrophotometer  
Model UV-1800  
Manufacturer Shimadzu  
Serial No. A11454908533 CD  
ID No. BKK\_EN0018  
Date of receipt 13 September 2024  
Date of calibration 13 September 2024  
Date of issue 13 SEP 2024

REVIEW BY	Juda K
APPROVED BY	Suk P
NEXT CAL DATE	13/09/2025

Customer name ALS Laboratory Group (Thailand) Co., Ltd.  
Address 104 Soi Phattanan 40, Phattanan Road, Phattanan, Suan Luang, Bangkok 10250  
Temperature (25.3 - 26.7) °C (On site)  
Humidity (50.4 - 55.9) %RH (On site)  
Equipment condition Good Operation  
Calibration Location Organic Preparation Lab  
Calibration Procedure In-house method WI-UV-702-01 based on ASTM E275-01  
Traceability Wavelength Accuracy is traceable to certificate No. 106372 and 106371  
Photometric Accuracy is traceable to certificate No. 106364 and 111398  
Stray Light is traceable to certificate No. 106377  
The above certificate are traceable to SI unit through Stama Scientific Ltd.  
(UKAS accredited calibration laboratory NO. 0659)  
Calibrated by Mr.Wanchana Janloy

Approved by

Mr.Sonthi Temboonsakdi  
Service Manager

The above results are valid exclusively for the calibrated item(s) as mention in this report / certificate.  
Advertising the report / Certificate and publicity of the results are prohibited and also shall not be reproduced  
except in full, without written approval of the Bara Scientific Co., Ltd.

FM-L13 108-30-05-57



Bara Scientific Co., Ltd.  
968 U Chu Liang Building Floor7 Ramat Road  
Silom Bangkok Bangkok Thailand 10500  
Tel : 02-6324300 Fax : 02-6375496-7  
www.barascientific.com



## Certificate of Calibration

Certificate No. BSCC-UV-374/24

Number of Page(s) 2 of 3

#### Calibration Results:

##### 1.Wavelength Accuracy

Certified Wavelength (nm)	UUC (nm)	Error (nm)	Uncertainty (±nm)
241.70	241.55	-0.15	0.18
334.02	333.85	-0.17	0.18
418.53	418.57	0.04	0.18
572.99	572.97	-0.02	0.18
879.41	879.17	-0.24	0.18

##### 2.Photometric Accuracy (UV)

Wavelength (nm)	Certified Absorbance (A)	UUC (A)	Error (A)	Uncertainty (±A)
235	0.0000	0.0000	0.0000	0.0075
	0.7171	0.7169	-0.0002	0.0075
257	0.0000	0.0000	0.0000	0.0075
	0.8354	0.8345	-0.0009	0.0075
313	0.0000	0.0000	0.0000	0.0075
	0.2786	0.2781	-0.0005	0.0075
350	0.0000	0.0000	0.0000	0.0075
	0.6199	0.6194	-0.0005	0.0075

\*CNR = Customer not request

The above results are valid exclusively for the calibrated item(s) as mention in this report / certificate.  
Advertising the report / Certificate and publicity of the results are prohibited and also shall not be reproduced  
except in full, without written approval of the Bara Scientific Co., Ltd.

FM-UV-708-02 Rev.01 (2301/63)



Bara Scientific Co., Ltd.  
968 U Chu Liang Building Floor7 Ramat Road  
Silom Bangkok Bangkok Thailand 10500  
Tel : 02-6324300 Fax : 02-6375496-7  
www.barascientific.com



## Certificate of Calibration

Certificate No. BSCC-UV-374/24

Number of Page(s) 3 of 3

#### Calibration Results:

##### 3.Photometric Accuracy (Visible)

Wavelength (nm)	Certified Absorbance (A)	UUC (A)	Error (A)	Uncertainty (±A)
420.0	0.0000	0.0000	0.0000	0.0042
	0.5781	0.5765	-0.0016	0.0042
	0.7119	0.7105	-0.0014	0.0042
	1.0189	1.0174	-0.0015	0.0042
440.0	0.0000	0.0000	0.0000	0.0042
	0.5610	0.5613	0.0003	0.0042
	0.7001	0.6984	-0.0017	0.0042
	1.0026	1.0011	-0.0015	0.0042
465.0	0.0000	0.0000	0.0000	0.0042
	0.5235	0.5232	-0.0003	0.0042
	0.6614	0.6598	-0.0016	0.0042
	0.9456	0.9444	-0.0012	0.0042
546.1	0.0000	0.0000	0.0000	0.0042
	0.5249	0.5245	-0.0004	0.0042
	0.6975	0.6956	-0.0019	0.0042
	1.0009	0.9994	-0.0015	0.0042
590.0	0.0000	0.0000	0.0000	0.0042
	0.5590	0.5586	-0.0004	0.0042
	0.7726	0.7708	-0.0017	0.0042
	1.1125	1.1114	-0.0011	0.0042
635.0	0.0000	0.0000	0.0000	0.0042
	0.5666	0.5666	0.0000	0.0042
	0.7620	0.7604	-0.0016	0.0042
	1.0962	1.0971	-0.0011	0.0042

\*CNR = Customer not request

##### 4.Stray Light

Standard cut-off wavelength (nm)	Wavelength (nm)	Transmission (%T)	Absorbance (A)
200.85±0.11nm	199.58	0.9520	2.0217

The Stray light transmission reference is less than 1.0%T and Stray light absorbance reference is greater than 2.00A  
\*Stray Light not NRC-CNRC Accredited.

The measurement uncertainty is base on a standard uncertainty multiplied by a coverage factor k=2 providing a level of confidence of approximately 95%.

\*\*\*End of Certificate\*\*\*

The above results are valid exclusively for the calibrated item(s) as mention in this report / certificate.  
Advertising the report / Certificate and publicity of the results are prohibited and also shall not be reproduced  
except in full, without written approval of the Bara Scientific Co., Ltd.

FM-UV-708-02 Rev.01 (2301/63)







Cert. No.: 23T2290  
Page: 2 of 2

#### Result of Calibration:-

Without Adjustment

Function: Temperature measurement  
Dimension of probe: Diameter 3 mm., Length 55 mm. Sheath material : Stainless Steel

Immersion Depth (mm.)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty of Measurement (±°C)
50	25.0003	24.9	-0.1003	0.12
50	30.0031	29.9	-0.1031	0.12
50	35.0003	35.0	-0.0003	0.12
50	40.0014	40.0	-0.0014	0.12
50	45.0029	45.0	-0.0029	0.12

UUC\* : Unit Under Calibration

The reported uncertainty of measurement was based on standard uncertainty multiplied by a coverage factor  $k = 2$ , providing a level of confidence of approximately 95%.

-000-

a 1195879



Agilent Technologies (Thailand) Limited  
U CHU LIANG BLDG. 22/F UNIT A.D.  
988 RAMA 4 ROAD, SILOMA BANGKOK  
Bangkok 10500 Thailand

Tel: +662 637 6363  
Fax: +662 637 6334  
Email: ccc-smt@agilent.com  
Website: www.agilent.com/cham

#### Customer Contact:


ALS Laboratory Group (Thailand) Co.  
Ltd Head Office  
104 Phatthanakan 40 Phatthanakan Rd  
Khuwaeng Phatthanakan Khet Suan  
TAX ID : 0105540004859  
chanattagarnimchom@alsglobal.com  
227158760

#### Invoice To:

ALS Laboratory Group (Thailand) Co.  
Ltd Head Office  
104 Phatthanakan 40 Phatthanakan Rd  
Khuwaeng Phatthanakan Khet Suan

#### SERVICE REPORT

Customer Purchase Order Number:	Customer Number: 70371013
Service Request:	Service Request Date:
Service Order: 6006676060	Service Confirmation: 6905905441

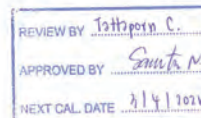

#### Delivery Site:

ALS Laboratory Group (Thailand) Co.  
Ltd Head Office  
104 Phatthanakan 40 Phatthanakan Rd  
Khuwaeng Phatthanakan Khet Suan

#### Location:

Room  
Bldg  
Lab  
Dept

#### Direct Inquiries to:

Contact Name: Customer Contact Center  
Contact E-mail: ccc-smt@agilent.com  
Contact Telephone: +662 637 6363  
Contact Fax: +662 632 4334

product applications software services

Learn more about Agilent's Special Offers, Products, Services and our full range of laboratory productivity solutions optimized for your applications and workflows. Visit us at [www.agilent.com/cham](http://www.agilent.com/cham)

Agilent Technologies (Thailand) Limited, Head Office  
U Chu Liang Bldg. 22/F Unit A.D.  
988 Rama 4 Road, Silom, Bangkok.  
Bangkok 10500 Thailand  
Tax ID : 0105542080218

DBank N.A. Bangkok Branch  
395 Interchange 21 Building, Sukhumvit Road, Klongroey New  
Sub-district, Wattana District, Bangkok 10110 Thailand  
Acc. No: 012-4402 007  
THB-Krung Thai Bank PCL  
Siem Square Bv. 416/1-2 Rama 1 Rd, Pathumwan, BKK 10330  
Thailand

Page 1 of 3

Service Confirmation Number: 6905905441  
Service Confirmation Date: 06.10.2024

Service Confirmation Number: 6905905441  
Service Confirmation Date: 06.10.2024

#### Service Instrument:

Model Number	Model Description	Serial Number	System Handle	Parent Asset
SYS-IM-7900	ICPMS 7900 System			
G8410A	SPS 4 Autosampler	AU15430722	ICP MS 7900	SYS-IM-7900
G8411A	ISIS 3 for Agilent 7850/7900/8900	JP15510227	ICP MS 7900	SYS-IM-7900
G3292A	PSC 8106T Chiller	ZU15A1948	ICP MS 7900	SYS-IM-7900
G8403A	Agilent 7900 ICP-MS	JP15471169	ICP MS 7900	SYS-IM-7900

#### Service Items:

Item	Service/Part #	Description	Qty	Entitlement	Service Start	Service End
1008	EQO	Enterprise Operational Qualification	1.00	Agreement Entitlement: 100 % covered	04.10.2024	04.10.2024
1010	6185-5950	ICP-MS Checkout Solutions	1.00	Agreement Entitlement: 100 % covered		

#### Additional Information:

#### Service Information:

<b>Problem Description:</b> *WU-EQO-IM-7900-6001253055		
<b>Service Provided:</b> Perform: OQ Hardware, Test CDS login, auto sampler, Auto tune, BG and 20 Min stability, I calibrate the instrument No BKK_EL0043 test all pass.		
<b>Service Overview Code:</b> Reason Code: Scheduled Service Diagnosis Code: Scheduled Service Resolution Code: Scheduled Service		
<b>Reported Hours:</b> 7.8	<b>Travel Hours:</b> 2.0	
<b>Customer Field Service Representative Name:</b> Pantkep Kuranathain	<b>Customer Field Service Representative Signature:</b> 	<b>Date:</b> 09 Oct 2024
<b>Customer Name:</b> Supakwan Mak	<b>Customer Signature:</b> 	<b>Date:</b> 08 Oct 2024
<b>Additional Comments:</b>		




Certificate No. T231676

Page 1 of 6

### Certificate of Calibration

Equipment : HEATING BLOCK  
Manufacturer : Environmental Express  
Model : SC 196  
Serial No. : 6974CECW3285  
Customer Code : BKK\_EL0054  
ID No. : T5306A3  
Customer : ALS Laboratory Group (Thailand) Co.,Ltd.  
104 Phatthanakan 40, Phatthanakan Rd., Khwaeng Phatthanakan,  
Khet Suan Luang, Bangkok 10250  
Customer Location : Acid Digestion Lab  
Date of Receipt : 13 September 2023  
Calibrated By : Saneek Musikanwan ( Site Calibration Manager )  
Approved By : / Sujjar Nakhakred ( Site Calibration Manager )  
Date of Issue : 26 SEP 2023

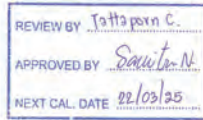

The uncertainties are for a confidence probability of approximately 95%.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standard laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the Metrological Center.

FM-L13 108/30-05-57



Certificate No. T231676

Page 2 of 6

### Calibration Report

Equipment : HEATING BLOCK  
Date of Calibration : 22 September 2023  
Environment : Temperature : 21.8-23.1 °C  
Line Voltage : 221.6-226.3 V  
Relative Humidity : 55 - 65 %RH

#### Condition of this results of calibration :

- This equipment was calibrated by insert 20 standard thermocouples type T into its chamber , the other one standard thermocouples type T use for ambient temperature measurement . The calibration was done in according to WI-T20.  
All data show below were final values and the initial data from customer request . The temperature scale used was based on ITS - 90 .
- Reference Standard Instrument :

Instrument	Model	Instrument No.	Certificate No.	Due Date
TC	TYPE T	TN21-TN30	T230014	17 January 2024
TC	TYPE T	TN31-TN40	T230014	17 January 2024
DATA LOGGER	34970A	T151	T230014	17 January 2024
- This certificate is traceable to :  
National Institute of Metrology ( Thailand ) through Metrological Center ( NSC-TISI-TIS 17025 CALIBRATION 0244.)
- Condition of calibrated item : good  
Equipment Description :  
Time Constant : 2 Hour 20 Minute At 95 °C  
Fresh Air Dumper : ☐ Open ☐ Min ☐ Medium ☐ Max  
☐ Close  
☒ Not Available
- Adjustment :  
( ) without adjustment ( X ) after adjustment

Approved By


FM-L13 108/30-05-57



Certificate No. T231676

Page 3 of 6

### Calibration Report

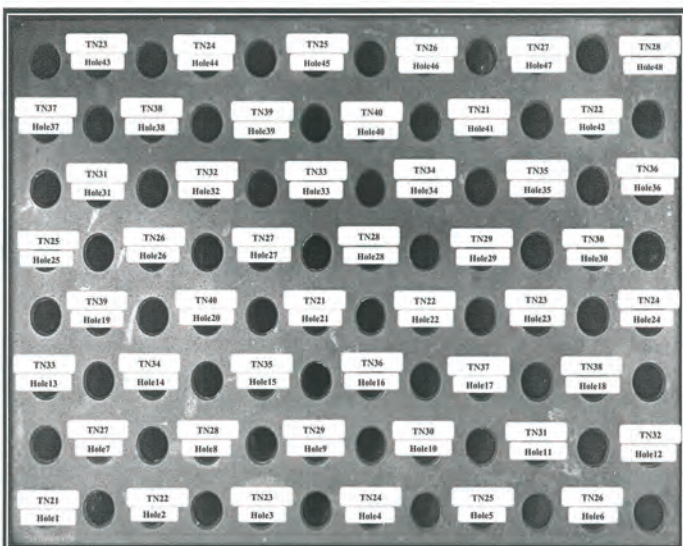

FRONT CONTROL

Approved By

FM-L13 108/30-05-57



Certificate No T231676

Page 4 of 6

### Calibration Report

#### Measurement Results

Calibration Point	Average Standard Reading at each position ( °C )					
R1 Hole1-Hole6	TN21	TN22	TN23	TN24	TN25	TN26
CAL POINT	Max	95.01	94.41	95.20	95.41	94.51
	Min	94.57	93.95	94.75	94.92	94.00
	Average	94.79	94.18	94.98	95.17	94.26
R2 Hole7-Hole12	TN27	TN28	TN29	TN30	TN31	TN32
	Max	95.36	95.43	95.19	95.16	95.35
	Min	94.94	94.95	94.72	94.71	94.90
	Average	95.15	95.19	94.96	94.94	95.13
R3 Hole13-Hole18	TN33	TN34	TN35	TN36	TN37	TN38
	Max	95.37	95.50	95.22	95.21	95.33
	Min	94.99	95.09	94.78	94.82	94.88
	Average	95.18	95.30	95.00	95.02	95.11
R4 Hole19-Hole24	TN39	TN40	TN21	TN22	TN23	TN24
	Max	95.59	94.42	94.52	94.24	94.63
	Min	95.21	94.06	94.13	93.88	94.28
	Average	95.40	94.24	94.33	94.06	94.45
R5 Hole25-Hole30	TN25	TN26	TN27	TN28	TN29	TN30
	Max	95.19	95.38	92.93	95.30	95.14
	Min	94.83	95.03	92.56	94.95	94.79
	Average	95.01	95.20	92.75	95.12	94.96
R6 Hole31-Hole36	TN31	TN32	TN33	TN34	TN35	TN36
	Max	94.63	94.90	94.77	94.31	94.24
	Min	94.24	94.55	94.44	93.98	93.92
	Average	94.43	94.72	94.60	94.14	94.08
R7 Hole37-Hole42	TN37	TN38	TN39	TN40	TN21	TN22
	Max	94.30	94.44	94.04	93.81	94.89
	Min	93.95	94.05	93.67	93.48	94.39
	Average	94.13	94.24	93.86	93.65	94.64
R8 Hole43-Hole48	TN23	TN24	TN25	TN26	TN27	TN28
	Max	95.99	95.63	95.28	95.29	95.45
	Min	95.57	95.13	94.82	94.84	94.99
	Average	95.78	95.39	95.05	95.07	95.22

Approved By

FM-L13 108/30-05-57





Certificate No T231676

Page 5 of 6

## Calibration Report

### Measurement Results

Calibration Point	Average Standard Reading at each position (°C)					
R1 Hole1-Hole6	TN21	TN22	TN23	TN24	TN25	TN26
CAL POINT	Max	105.23	104.32	105.43	105.25	104.44
105	Min	104.94	105.95	105.15	105.04	104.11
	Average	105.09	104.13	105.29	105.15	104.28
R2 Hole7-Hole12	TN27	TN28	TN29	TN30	TN31	TN32
	Max	105.30	105.12	105.18	105.22	105.12
	Min	105.11	104.92	104.96	105.00	104.92
	Average	105.20	105.02	105.07	105.11	105.02
R3 Hole13-Hole18	TN33	TN34	TN35	TN36	TN37	TN38
	Max	105.37	105.63	105.02	104.80	104.69
	Min	105.17	105.37	104.75	104.59	104.50
	Average	105.27	105.50	104.88	104.69	104.60
R4 Hole19-Hole24	TN39	TN40	TN21	TN22	TN23	TN24
	Max	105.31	104.43	106.41	104.71	105.63
	Min	105.08	104.22	106.15	104.81	105.37
	Average	105.19	104.33	106.28	104.56	105.50
R5 Hole25-Hole30	TN25	TN26	TN27	TN28	TN29	TN30
	Max	104.95	106.26	103.34	105.78	105.59
	Min	104.67	105.96	103.08	105.56	105.36
	Average	104.81	106.11	103.21	105.67	105.48
R6 Hole31-Hole36	TN31	TN32	TN33	TN34	TN35	TN36
	Max	104.75	104.86	104.80	105.20	104.50
	Min	104.54	104.63	104.59	105.00	104.32
	Average	104.65	104.75	104.69	105.10	104.41
R7 Hole37-Hole42	TN37	TN38	TN39	TN40	TN21	TN22
	Max	104.30	104.90	104.85	104.65	104.88
	Min	104.09	104.72	104.66	104.49	104.63
	Average	104.19	104.81	104.75	104.57	104.76
R8 Hole43-Hole48	TN23	TN24	TN25	TN26	TN27	TN28
	Max	105.71	105.85	105.39	105.61	105.42
	Min	105.45	105.61	105.14	105.27	105.18
	Average	105.58	105.73	105.27	105.44	105.30

Approved By. \_\_\_\_\_

FM-L13 108/30-05-57



Certificate No. T231676

Page 6 of 6

## Calibration Report

### Measurement Results:

HEATING BLOCK			Temperature Distribution	
Setting (°C)	Reading (°C)		Stability (±°C)	Uncertainty (±°C)
	Min, Max	Average		
100.0	100.3, 100.5	100.4	0.26	0.81
107.0	107.0, 107.1	107.1	0.19	0.78

\* The quoted uncertainty exclude " uniformity "

The calibration result apply only the above calibrated item.

The result of test was found accurate as shown on date and place of test only.

The reported expanded uncertainty is based on a standard uncertainty multiplied by a coverage factor k=2 which for a t-distribution, providing a level of confidence of approximately 95 %.

Approved By. \_\_\_\_\_

FM-L13 108/30-05-57

BKK\_EL0037



Agilent CrossLab Start Up Services

## Agilent 5100 5110 ICP-OES Preventive Maintenance

REVIEW BY	Thirina B.
APPROVED BY	Sana L.N.
NEXT CAL DATE	23/02/2025

Agilent Preventive Maintenance provides factory recommended service for your analytical instruments to assure reliable operation and the accuracy of your results

Delivered by highly trained and certified service engineers using genuine Agilent parts and supplies, Agilent Preventive Maintenance provides what you need to reduce unplanned downtime and keep your systems operating at their peak performance.

This checklist is used as a guide for completing the preventive maintenance tasks. A signed copy of this checklist is provided for your records.



Agilent 5100, 5110 Preventive Maintenance Checklist



## Introduction

### Customer Information

- Customers should provide all necessary operating supplies upon request of the engineer.
- A customer representative should be available to the engineer while performing the preventive maintenance procedures. Customers are responsible for regular maintenance and are encouraged to observe the service representative.
- Any parts not included in the Parts Lists section of this document are not part of the recommended Preventive Maintenance service nor are they included in the price of this service.
- If a system requires the use of extra or special procedures and/or parts for the maintenance service, then these must be ordered separately and charged as a repair, which may incur additional costs.
- For customers using HIF applications, the instrument should be returned to its standard sample introduction system.



### Important Customer Web Links

- To access **Agilent University**, visit <http://www.agilent.com/crosslab/university/> to learn about training options, which include online, classroom and onsite delivery. A training specialist can work directly with you to help determine your best options.
- To access the **Agilent Resource Center** web page, visit <https://www.agilent.com/en-us/agilentresources>. The following information topics are available:
  - Sample Prep and Containment
  - Chemical Standards
  - Analysis
  - Service and Support
  - Application Workflows
- The **Agilent Community** is an excellent place to get answers, collaborate with others about applications and Agilent products, and find in-depth documents and videos relevant to Agilent technologies. Visit <https://community.agilent.com/welcome>
- Videos about specific preparation requirements for your instrument can be found by searching the **Agilent YouTube** channel at <https://www.youtube.com/user/agilent>
- Need to place a service call?** Flexible Repair Options | Agilent

### Service Engineer's Responsibilities

- Contact the customer and ensure that all necessary supplies are available before the preventive maintenance visit.
- Only select those pages that relate to the system or module being serviced.
- Complete empty fields with the relevant information.
- Complete the relevant checkboxes in the checklist using either a "X" or tick mark "✓".
- Check "Service not applicable" check boxes to indicate services/tasks not delivered, as appropriate.
- Complete the Preventive Maintenance services in the most logical order relevant to the individual system service in the order of the tasks listed.
- Complete the **Service Review** section together with the customer.
- Complete the fields for page numbers at the foot of each selected page.
- Add relevant page numbers to selected pages and complete the total number of pages field in the Service Completion section.
- Ask the customer to sign the Service Verification section including the customer's and your signature.**

## Instrument Maintenance

### System Information

- ☒ Check this box if an instrument configuration report is attached instead of completing the table:

Instrument System Name and ID	2401A / MY1610005
Instrument System Site and Location	RLS Laboratory Group (Thailand) Co., LTD.

List System Component Product Numbers	List the Serial Numbers of each Component
1. 873010A	MY16010005
2. 84101A	AU1544-0764
3. 82411-85001	1004-00169
4.	
5.	
6.	
7.	
8.	
9.	

ICP-OES Configuration Table	Circle the type or write in the type if other
Nebulizer Type	Spray / On/Off / Conical / Other
Spray Chamber	Cyclonic Single Pass / Cyclonic Double Pass / Other
Torch	Radial / Dual View / Other
Torch Type	One Piece / Semi Detachable / Fully Detachable / Other
Injector Diameter	2.4mm / 1.8mm / 1.4mm / 0.8mm / Other
Injector Material	Quartz / Ceramic / Other

### Preparation

- ☒ Discuss any specific issues with the customer before starting.
- ☒ Review the instrument logbook for recorded problems and comments.
- ☒ Save instrument control settings before starting the procedure.
- ☒ Perform a general inspection of the system for cleanliness.
- ☒ Check for proper installation of parts, assemblies, sensors etc.
- ☒ Check system for required installation of components and implementation of Service Notes.
- ☒ Check for required firmware/software updates and verify with customers if they would like them installed.
- ☒ For HF application systems, if standard sample introduction system was not installed, ask the customer to install it.
- ☒ Ask the customer to remove any samples from the ICP-OES sample introduction area, auto sampler or around the ICP-OES.



## Preventive Maintenance Procedures

### Record Pre-PM instrument performance

- ☒ Run Instrument Performance test
- ☒ Record results in Instrument Performance Test Results Table – Pre-PM.

### Clean and inspect ICP-OES system

- ☒ Look for any obvious external damage or problems.
- ☒ Inspect water cooling hoses, gas lines and power cord for excessive wear or damage.
- ☒ Perform a general internal inspection of the system for excessive dust accumulation, clean if necessary.
- ☒ Inspect sample introduction components and record any required maintenance in the Service Engineer Comments and notify the customer as the required actions required.
- ☒ Record the instrument operating conditions in the ICP-OES Status Results Table.
- ☒ Replace the polychromator purge filter.
- ☒ Replace the radial pre-optics window
- ☒ Replace the axial pre-optics window for SVDV and VDV instruments.
- ☒ Check exhaust flow for the correct positive extraction at the exhaust duct to insure they meet minimum specifications.
- ☒ Replace air inlet dust filter.
- ☒ Replace high capacity air inlet dust filter element if installed.
- ☒ Remove and clean instrument water inlet filter.

### Agilent Water Recirculator

- ☐ Service not applicable
- ☒ Drain cooling fluid and remove any particles from the chiller reservoir
- ☒ Remove, clean and reinstall water inlet metal mesh filter if present
- ☒ Re fill with Agilent Cool Clear cooling fluid
- ☒ Clean the cooling system Air filter and the condenser.

### SPS 3 Auto Sampler

- ☒ Service not applicable
- ☐ Power cycle the autosampler and verify successful initialization.
- ☐ Inspect X and Z axis belts for wear. Replace is necessary.
- ☐ Clean X and Z axis slide shafts.
- ☐ Using customer's racks and the Agilent software move the sample probe to the 4 outermost corners and rinse port, ensure that the probe is approximately centered in the vial.

### SPS 4 Auto sampler

- ☐ Service not applicable
- ☒ Clean the spill tray, rack location mat, end frames and chassis with a damp soft cloth and diluted mild detergent.
- ☒ Clean the auto sampler cover panels, if cover kit is installed, with domestic window cleaner.
- ☒ Check the X-axis and Z-axis drive belts for cracks, splits, damaged teeth, excessive fraying, color changes or degradation from fumes.
- ☒ Check the X-axis, Theta-axis and Z-axis FFC cables for cracks, incorrect positioning, damaged edges or damaged connectors.
- ☒ Pump Tubing Replacement: Replace peristaltic pump tubing. Replace all tubing that goes from the rinse station to the pump and from the pump to the waste/rinse bottles
- ☒ Test using customer's tray and move the sample probe to the sample vial 1, wash vial and rinse port and ensure that the probe is centered in the vial. If not use calibration wizard and calibrate the position.

### AVS 4, 6, 7 Advanced Valve System

- ☒ Service not applicable
- ☐ Replace valve rotor seal
- ☐ Check fittings for signs of leaks
- ☐ Check tubing including autosampler tubing for kinks or excessive wear
- ☐ Check high flow pump for signs of leaks

### ICP-OES adjustment

- ☒ Check position of Zn peak, adjust if required.
- ☒ Check Argon Ratio, adjust to specified value if required.
- ☒ Perform Detector Calibration.
- ☒ Perform Instrument Calibration

### Record Post-PM instrument performance

- ☒ Run Instrument Performance test.
- ☒ Record results in Instrument Performance Test Results Table - Post PM
- ☒ For systems using ICP Expert version 7.3 and above, run the following Instrument tests
  - ☒ Subsystem Communications Test
  - ☒ Air Flow
  - ☒ Water Flow
  - ☒ Gas Flows
  - ☒ RF Generator
  - ☒ Camera Test
  - ☒ Optics Test
  - ☒ Nebulizer Test
- ☒ Record the result in the Instrument Test Results Table

### Restore Instrument

- ☒ For HF applications, ask the customer to reinstall their sample introduction system.
- ☒ Leave system in an idle state: on and purging.
- ☒ Guidance: If the PM service is performed prior to a qualification service, then use the qualification procedure as a guide for final instrument set up and checkout.

### Service Review

- ☒ Attach available reports/printouts of all tests to this documentation.
- ☒ Record the Preventive Maintenance service activity in the customer's records/logbook.
- ☒ Record the PM event in the Smart Alerts logbook, if applicable.
- ☒ Update/reset instrument maintenance counters as appropriate
- ☒ Affix the PM sticker to the system or instrument logbook based on the customer's request.
- ☒ Complete the Service Engineer Comments section if there are additional comments.
- ☒ Review this service, parts replaced, and test results obtained with the customer.
- ☒ If the instrument firmware was updated, record the details of the change in the Service Engineer's Comments box. Systems in a compliant environment may need additional documentation.
- ☒ Complete the Signature Page with both Service Engineer and Customer signatures.

## Test Results

### Instrument Performance Test Results Table

Note: These measurements do not form part of any specification and are for reference only.

	Pre PM Sensitivity Check		Post PM Sensitivity Check	
	Radial	Axial*	Radial	Axial*
Zn 213.857 nm SBR	1511.1	2444.3	1520.0	2421.8
Mn 257.610 nm SBR	1351.1	1956.6	2348.5	1999.3
Al 396.152 nm SBR	1.1	116.0	5.8	10.3
K 766.491 nm SBR	3.5	64.0	5.6	92.2

\* Axial result is not applicable for G8016AA, G8012AA Radial View Instruments

### Instrument Test Results Table

Note: The Instrument Test results are for systems using ICP Expert version 7.3 and above only.

Instrument Test	Result
Subsystem Communications Test	PASS
Air Flow	PASS
Water Flow	PASS
Gas Flows	PASS
RF Generator	PASS
Camera Test	PASS
Optics Test	PASS
Nebulizer Test	PASS

## ICP-OES Status Results Table

Note: These measurements do not form part of any specification and are for reference only.

Measurement	Standby Mode	Plasma On
Main Voltage	219.31 V	VAC 215.12 V
Main Current	0.09 L	0.113 A
Instrument Temperature	11.8 °C	23.8 °C
RF Air Flow (sensor speed)	94.0 Hz	23.0 Hz
Plasma Exhaust Temperature	No measurement	50.1 °C
Water Flow Oscillator	No measurement	1.20 L/min
Water Flow Detector	1.14 L/min	1.09 L/min
Water Inlet Temperature	27.8 °C	23.6 °C
Polychromator Temperature	35.0 °C	36.0 °C
CCD Temperature	-40.1 °C	-40.0 °C
Thermal Stabilizer	31.8 °C	34.6 °C
Argon Supply Pressure	514.94 kPa	557.10 kPa
Purge Gas Supply Pressure*1	610.61 kPa	514.30 kPa
Option Gas Supply Pressure*1	— kPa	— kPa
Nebulizer Flow	No measurement	0.70 L/min
Nebulizer Back Pressure	No measurement	236.06 kPa
Plasma Gas Flow	No measurement	11.89 L/min
Auxiliary Gas Flow	No measurement	1.00 L/min
RF Power	No measurement	1599.6 W
RF Supply Current	No measurement	5.66 A
RF Supply Voltage	No measurement	164.81 V

\*1 If option installed

## Consumed PM Parts

Part Description	Part Number	Product or Model# where used	Quantity consumed
Axial Pre-Optic Window	G8010-68014	G8010A, G8011A, G8014A/G8015A	1
Radial Pre-Optic Window	G8010-68015	All	1
Agilent Cool Clear Coolant Fluid	5799-0037	Agilent Water Recirculator	1
Purge Gas Filter	G8010-60136	All	1
Air inlet filter	G8000-68002	All	1
High Capacity Air Filter	G8010-60189	Optional	—
Rotor seal for 6-port valve for AYS6/7	G8494-60002	G8494A/G8495	—
Rotor seal for 4-port valve for AYS4	G8492-60002	G8493A	—
Rinse solution to rinse station 2.5mm ID x 1m	G8410-80123	SPS 4	1
Barb connector 2.5mm x 1.5mm ID	G8410-80124	SPS 4	1
PVC waste tubing 1/8mm OD x 5mm ID, 2m	G8410-80122	SPS 4	1

### Additional Parts may be required from engineer's stock:

X axis drive belt	5410047800	SPS 3	—
Z axis drive belt	5410047400	SPS 3	—
Peristaltic pump tubing (PVC SoluFlex 3 bridge)	3710049000	SPS 4	—

### Consumed Parts Reference

(Purchased by customer, not included as part of PM)

Section Not Applicable

Part Description	Part Number	Product or Model# where used	Quantity consumed
------------------	-------------	------------------------------	-------------------

## Signature Page


### Service Engineer Comments (optional)

If there are any specific points you wish to note as part of performing the installation or other items of interest for the customer, please write in this box.

### Service Verification

Service Request Number: 6006602534  
Date Service Completed: Feb 21, 2022  
Service Engineer Name: Nelson, Lucretia  
Customer Name:   
Service Engineer Signature:   
Customer Signature:   
Total number of pages in this document: 14



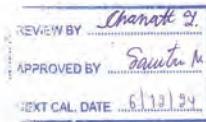

### Performance Verification Certificate for Mercury Analyzer

**PRODUCT ID** Quicktrace M-8000 , Teledyne Leeman Labs  
**Equipment ID** BKK\_EL0128 Mercury Analyzer  
 S/N: US22133002  
 BKK\_EL0129 Autosampler  
 S/N: 052222A560  
**Customer Name** ALS Laboratory Group (Thailand) Co., Ltd.  
**Address** 104 Soi Pattana 40, Pattana Rd. Suan Luang, Suan Luang  
 Bangkok 10250 Thailand

**Date of Qualified** December 6, 2023  
**Next Due date** December 6, 2024

This certifies for products which was performed in acceptable criteria specifications

<b>Autosampler &amp; Sample Introduction</b>	<b>PASSED</b>
<b>Analyzer</b>	<b>PASSED</b>
<b>Gas Liquid Separator &amp; Dryer</b>	<b>PASSED</b>
<b>CVAFS Detector</b>	<b>PASSED</b>
<b>Electronics/Mechanical</b>	<b>PASSED</b>
<b>Data station/PC</b>	<b>PASSED</b>
<b>Analytical test</b>	<b>PASSED</b>

Provided by

Scientist Instrument Co.,Ltd.  
 113 Soi Ekachai 44, Ekachai Road  
 Khlong Bang Phran, Bangbon  
 Bangkok 10150 Thailand

**Certified by**   
 Thunraphol Sakdayos  
 Service Engineer



### Performance Verification Certificate for Mercury Analyzer

**PRODUCT ID** Quicktrace M-8000 , Teledyne Leeman Labs  
**Equipment ID** BKK\_EL0128 Mercury Analyzer  
 S/N: US22133002  
 BKK\_EL0129 Autosampler  
 S/N: 052222A560  
**Customer Name** ALS Laboratory Group (Thailand) Co., Ltd.  
**Address** 104 Soi Pattana 40, Pattana Rd. Suan Luang, Suan Luang  
 Bangkok 10250 Thailand

**Date of Qualified** December 6, 2024  
**Next Due date** December 6, 2025

This certifies for products which was performed in acceptable criteria specifications

<b>Autosampler &amp; Sample Introduction</b>	<b>PASSED</b>
<b>Analyzer</b>	<b>PASSED</b>
<b>Gas Liquid Separator &amp; Dryer</b>	<b>PASSED</b>
<b>CVAFS Detector</b>	<b>PASSED</b>
<b>Electronics/Mechanical</b>	<b>PASSED</b>
<b>Data station/PC</b>	<b>PASSED</b>
<b>Analytical test</b>	<b>PASSED</b>

Provided by

Scientist Instrument Co.,Ltd.  
 113 Soi Ekachai 44, Ekachai Road  
 Khlong Bang Phran, Bangbon  
 Bangkok 10150 Thailand

**Certified by**   
 Thunraphol Sakdayos  
 Service Engineer




TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
 CORPORATE SERVICES & EQUIPMENT CALIBRATION AND TESTING SERVICES  
 534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250  
 TEL: 0-2717-3000-29 FAX: 0-2719-9484



**Cert. No.:** 23TM1103  
**Page:** 1 of 3

### Certificate of Calibration

**Equipment :** Autoclave  
**Manufacturer :** Sanyo  
**Model :** MLS-3781  
**Serial No. :** 830167  
**ID No. :** BKK\_ML0037

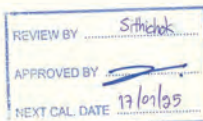

**Submitted by :** ALS Laboratory Group (Thailand) Co.,Ltd.  
 104 Phatthanakan 40, Phatthanakan Rd.,  
 Khwaeng Phatthanakan, Khet Suan Luang,  
 Bangkok 10250 Thailand

**Location :** Media Preparation Room

**Received Order :** 17 July 2023  
**Calibration Date :** 17 July 2023  
**Ambient Temperature :** ( 26 ± 10 ) °C  
**Relative Humidity :** ( 50 ± 30 ) %

**Calibrated by :** Preecha Hlahib

**Approved by :**   
 Approved Signatory

( ) Pornthippa Tameyakul  
 ( ) Malee Butkruea  
 (✓) Suwit Imjai

**Issue Date :** 24 July 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
 Approval of the head of Corporate Services : Equipment Calibration and Testing Services.



**Equipment :** Autoclave  
**Condition As-Received :** Used Item  
**Reference :** 2307-0285OC-3  
**Procedure Used :-**

**Cert. No.:** 23TM1103  
**Page:** 2 of 3

Calibration were conducted using in-house calibration procedure CP-OT03 according to direct measurement method with Data Acquisition which connected with Thermocouple Type T  
 The temperature scale used was based on ITS-90.

#### Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1 ) Data Acquisition	MY57013823	23LM66	TPA	25 Mar 2024

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

4. This result of calibration covers laboratory autoclaves for the sterilization of goods and material which

could be infected with organisms categorized as Hazard Group 1, 2 and 3\*\*

(\*\* = Categorization of pathogens according to hazard and categories of containment, second edition, 1990 )

It does not cover autoclaves for use with material infect with organisms in Hazard Group 4, for which complete containment and sterilization of infected condensate is considered to be essential.

This result of calibration does not apply to sterilizers or disinfectors used for medical, dental, pharmaceutical or veterinary purposes which are directly concerned with patient care, or those used for fabrics subjected to sterilization which are required to be dry at the end of cycle.

**Remark :** TPA : Technology Promotion Association ( Thailand - Japan )

**Result of Calibration :-** ( \* ) Without Adjustment

**Function of UUC\* :** Temperature Source



	Environmental		
	( °C )	( %R.H. )	( Volt )
Beginning of Calibration	22	53	220
Finished of Calibration	22	54	220

Position	Description	Ref. Std. ID No.:
1 =	Center of chamber	22-17TC-01
2 =	Temperature sensor	23-17TC-02
3 =	Exhaust port	19-17TC-03





Equipment : Autoclave  
Condition As-Received : Used Item  
Reference : 2307-0285OC-3  
Result of Calibration :- ( \* ) Without Adjustment  
Function of UUC\* : Temperature Source

Cert. No.: 23TM1103  
Page : 3 of 3

Operating parameter Set : Temperature = 121 °C  
Sterilization period = 15 minute

UUC* Setting ( °C )	UUC* Reading ( °C )	Position	Average* Standard Reading ( °C )	Stability ( ± °C )	Pressure Reading ( MPa )	Uncertainty ( ± °C )	Coverage Factor k
121	121	1	120.877	0.39	0.12	1.0	2
		2	120.870				
		3	120.866				

Average\* : The average of 30 values in each position.

Stability : One-half of the greatest maximum difference of measured temperature at any one probe.

UUC\* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-000-

a 1159504



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)  
CORPORATE SERVICES & EQUIPMENT CALIBRATION AND TESTING SERVICES  
534/4 PATTANAKARN ROAD 50/1 H. SUANLUANG, SUANLUANG BANGKOK 10250  
TEL: 0-2717-3000-29 FAX: 0-2719-9184



Cert. No.: 23TM1146  
Page : 1 of 3

## Certificate of Calibration

Equipment : Incubator

Manufacturer : SHEL-LAB

Model : 1915A

Serial No. : 0200599

ID No. : BKK\_ML0010

Submitted by : ALS Laboratory Group (Thailand) Co., Ltd.  
104 Phatthanakan 40, Phatthanakan Rd.,  
Khwaeng Phatthanakan, Khet Suan Luang,  
Bangkok 10250 Thailand

Location : Incubation & Micrological Reading

Received Order : 17 July 2023

Calibration Date : 17 July 2023

Ambient Temperature : ( 26 ± 10 ) °C

Relative Humidity : ( 50 ± 30 ) %

Calibrated by : Man Pattanapongpaiboon

Approved by :   
Approved Signatory

( ) Pongthipha Tameyaskul  
( / ) Malee Butkruea  
( ) Suwit Imjai

Issue Date : 24 July 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3 : Equipment, Calibration and Testing Services

A 0056489



Equipment : Incubator  
Condition As-Received : Used Item  
Reference : 2307-0285OC-1  
Procedure Used :-

Cert. No.: 23TM1146  
Page : 2 of 3

Calibration were conducted using calibration procedure CP-OT02 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector ( RTD ).

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1 ) Data Acquisition	MY49001451	23LM27	TPA	25 Feb 2024

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

Remark : TPA : Technology Promotion Association ( Thailand - Japan )


Result of Calibration :- ( \* ) Without Adjustment

Function of UUC\* : Temperature Source

Fresh air setting : Close

Environment during calibration		
	Beginning	Finished
Temp. ( °C )	24	24
REL Humid. ( % )	54	56
AC Supply ( Volt )	221	223

Position :	Ref. Std. ID No.:
1	19RTD-2/1
2	19RTD-2/2
3	19RTD-2/3
4	19RTD-2/4
5	19RTD-2/5
6	19RTD-2/6
7	19RTD-2/7
8	19RTD-2/8
9 (ref.)	19RTD-2/9

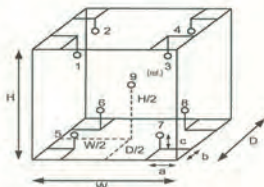

Probe Installation Details :

a = 10 cm  
b = 10 cm  
c = 10 cm

Dimension of Chamber :

D = 0.50 m  
W = 0.75 m  
H = 1.2 m  
Capacity = 0.45 m<sup>3</sup>

a 1172189



Equipment : Incubator  
Condition As-Received : Used Item  
Reference : 2307-0285OC-1  
Result of Calibration :- ( \* ) Without Adjustment  
Function of UUC\* : Temperature Source  
Fresh air setting : Close

Cert. No.: 23TM1146  
Page : 3 of 3

Calibration Point ( °C )	UUC* Setting ( °C )	UUC* Reading ( °C )	Temperature stability ( ± °C )	Temperature uniformity ( °C )	Overall Variation ( °C )	Coverage Factor k
35.0	35.0	35.0	0.055	0.30	0.44	2

Calibration Point ( °C )	Measured Temperature ( °C )									Uncertainty ( ± °C )
	1	2	3	4	5	6	7	8	9 (ref.)	
35.0	34.888	34.933	34.813	34.813	35.064	35.019	35.156	35.141	35.087	0.30

Average\* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC\* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-000-

a 1172188





## Certificate of Calibration

Cert. No.: 24TM667  
Page : 1 of 3

Equipment : Hot Air Oven  
Manufacturer : Binder  
Model : ED 240/E2  
Serial No. : 00-15533  
ID No. : BKK\_ML0013

Submitted by : ALS Laboratory Group (Thailand) Co.,Ltd.  
104 Phatthanakan 40, Phatthanakan Rd.,  
Khaeng Phatthanakan, Khet Suan Luang,  
Bangkok 10250 Thailand  
Location : Media Preparation Room

Received Order : 23 April 2024  
Calibration Date : 23 April 2024  
Ambient Temperature :  $(26 \pm 10) ^\circ\text{C}$   
Relative Humidity :  $(50 \pm 30) \%$

Calibrated by : Tawatchai Pama

Approved by :

( ) Ponpan Paipim  
(✓) Suwit Imjai  
( ) Kunchit Promprat

Issue Date : 26 April 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Equipment : Hot Air Oven  
Condition As-Received : Used Item  
Reference : 2404-0439OC-8  
Procedure Used :-

Cert. No.: 24TM667  
Page : 2 of 3

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Thermocouple Type T.  
The temperature scale used was based on ITS-90.

### Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1 ) Data Acquisition	MY49001451	24LM44	TPA	17 Mar 2025

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.


Remark : TPA : Technology Promotion Association ( Thailand - Japan )

Result of Calibration :- ( \* ) Without Adjustment

Function of UUC\* : Temperature Source

Fresh air setting : Close

Environment during calibration		
	Beginning	Finished
Temp. ( °C )	24	23
REL.Humid. ( % )	65	65
AC Supply ( Volt )	223	222

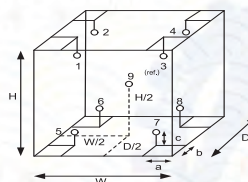

### Probe Installation Details :

a = 10 cm  
b = 10 cm  
c = 10 cm

### Dimension of Chamber :

D = 0.50 m  
W = 0.80 m  
H = 0.60 m  
Capacity = 0.24 m<sup>3</sup>

Position :	Ref. Std. ID No.:
1	24-19TC-01
2	24-19TC-02
3	24-19TC-03
4	24-19TC-04
5	24-19TC-05
6	24-19TC-06
7	24-19TC-07
8	24-19TC-08
9 (ref.)	24-19TC-09



Equipment : Hot Air Oven  
Condition As-Received : Used Item  
Reference : 2404-0439OC-8  
Result of Calibration :- ( \* ) Without Adjustment  
Function of UUC\* : Temperature Source  
Fresh air setting : Close

Cert. No.: 24TM667  
Page : 3 of 3

Calibration Point ( °C )	UUC* Setting ( °C )	UUC* Reading ( °C )	Temperature stability ( ± °C )	Temperature uniformity ( °C )	Overall Variation ( °C )	Coverage Factor k
180	180	180	0.64	2.7	3.7	2

Calibration Point ( °C )	Measured Temperature ( °C )									Uncertainty ( ± °C )
	1	2	3	4	5	6	7	8	9 (ref.)	
180	181.009	181.511	180.922	181.359	181.217	183.659	181.664	181.986	181.474	1.5

Average\* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC\* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-o-o-



## Certificate of Calibration

Cert. No.: 24TM469  
Page : 1 of 3

Equipment : Water Bath  
Manufacturer : Memmert  
Model : WNE 45  
Serial No. : L712.0429  
ID No. : BKK\_ML0056

Submitted by : ALS Laboratory Group (Thailand) Co.,Ltd.  
104 Phatthanakan 40, Phatthanakan Rd.,  
Khaeng Phatthanakan, Khet Suan Luang,  
Bangkok 10250 Thailand

Location : Incubation \$ Microbiological Reading

Received Order : 01 March 2024  
Calibration Date : 01 March 2024  
Ambient Temperature :  $(26 \pm 10) ^\circ\text{C}$   
Relative Humidity :  $(50 \pm 30) \%$

Calibrated by : Krisda Malee

Approved by :

( ) Ponthippa Tameyakul  
( ) Unnopphol Harachai  
(✓) Suwit Imjai

Issue Date : 4 March 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written  
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Equipment : Water Bath  
Condition As-Received : Used Item  
Reference : 2403-0001OC-1  
Cert. No.: 24TM469  
Page : 2 of 3

Procedure Used :-  
Calibration were conducted using in-house calibration procedure CP-OT04 Based on ASTM E715 according to direct measurement method with Data Acquisition which connected with Industrial Platinum Resistance Thermometer ( IPRT ).

The temperature scale used was based on ITS-90.

**Condition of this result of calibration**

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1 ) Data Acquisition	MY57013711	23LM115	TPA	11 Jul 2024

2. This certificate is valid only to the item calibrated on date and place of calibration.  
3. This certification is traceable to the International System of Unit.


Remark : TPA : Technology Promotion Association ( Thailand - Japan )

Result of Calibration :- ( \* ) Without Adjustment

Function of UUC\* : Temperature Source

Heat transfer medium used : Water

	Environmental		AC Voltage Supply
	( °C )	( %R.H. )	( Volt )
Beginning of Calibration	24	55	221
Finished of Calibration	23	56	220

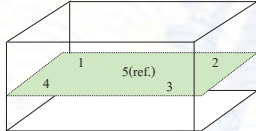

Position :	Ref. Std. ID No.:
1	4803988-001
2	4803988-002
3	4803988-003
4	4803988-004
5(ref.)	4803988-005



Equipment : Water Bath  
Condition As-Received : Used Item  
Reference : 2403-0001OC-1  
Cert. No.: 24TM469  
Page : 3 of 3  
Result of Calibration :- ( \* ) Without Adjustment  
Function of UUC\* : Temperature Source

Calibration point ( °C )	UUC* Setting ( °C )	UUC* Reading ( °C )	Average* Standard Reading ( °C )					Uncertainty  ( ± °C )
			Position					
			1	2	3	4	5 (ref.)	
44.5	44.5	44.5	44.469	44.462	44.492	44.510	44.496	0.15
45.0	45.0	45.0	44.975	44.974	45.007	45.023	44.999	0.15

Calibration point ( °C )	Uniformity ( °C )	Stability ( ± °C )	Coverage Factor k
44.5	0.087	0.029	2
45.0	0.069	0.031	2

Average\* : The average of 30 values in each position.

Uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Stability : One-half of the greatest maximum difference of measured temperature at any one probe.

UUC\* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity.

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-o0o-



## ภาคผนวก ด

สำเนาหนังสืออนุญาตขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

ที่ อก ๐๓๑๐(๑)/ ๑๖๑๖๘



กรมโรงงานอุตสาหกรรม  
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท  
เขตราชเทวี กรุงเทพฯ ๑๐๕๐๐

๒๐ พฤศจิกายน ๒๕๖๖

เรื่อง ค่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

เรียน กรรมการผู้จัดการ บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารเคมีของห้องปฏิบัติการวิเคราะห์เอกชน ลงวันที่ ๔ สิงหาคม ๒๕๖๖

สิ่งที่ส่งมาด้วย ๑. รายชื่อผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๑ แผน

๒. รายชื่อเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๕ แผน

๓. ขอบข่ายสารเคมีที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๓๑ แผน

ตามหนังสือที่อ้างถึง บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด ขอต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ๖-๒๐๔ สดงานที่ตั้งเลขที่ ๑๐๔ ซอยพัฒนาการ ๔๐ ถนนพัฒนาการ แขวงพัฒนาการ เขตสวนหลวง กรุงเทพมหานคร ต่อกรมโรงงานอุตสาหกรรม นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ไม่พบข้อบกพร่องของเอกสารประกอบคำขอต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน โดยมีองค์ประกอบดังนี้

ก. ผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๖ ราย ตามสิ่งที่ส่งมาด้วย ๑

ข. เจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ ๑๘๑ ราย ตามสิ่งที่ส่งมาด้วย ๒

ค. ขอบข่ายสารเคมีที่ได้รับขึ้นทะเบียนไม่วิเคราะห์ในน้ำเสีย น้ำใต้ดิน อากาศเสีย สิ่งปฏิกูล หรือวัตถุที่ไม่ใช่แล้ว และดิน ตามสิ่งที่ส่งมาด้วย ๓

หนังสือฉบับนี้จะหมดอายุในวันที่ ๒ กันยายน ๒๕๖๙ หากประสงค์จะต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ให้ยื่นคำขอต่ออายุพร้อมเอกสารประกอบคำขอต่อกกรมโรงงานอุตสาหกรรม ภายใน ๓๐ วัน ก่อนวันสิ้นสุดของหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ทั้งนี้ สามารถยื่นคำขอผ่านระบบอิเล็กทรอนิกส์ได้ที่หน้าเว็บไซต์กรมโรงงานอุตสาหกรรม

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

(นายศิระ จันทโรจน์)

นักวิทยาศาสตร์เชี่ยวชาญ วิชาการการแพทย์  
ผู้อำนวยการกองวิจัยและประเมินผลพิษโรงงาน  
ปฎิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม

กองวิจัยและประเมินผลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕

โทรสาร ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓

ไปรษณีย์อิเล็กทรอนิกส์ sarabang@diw.mail.go.th



"อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว"



เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด เลขทะเบียน ๖-๒๐๔

ที่ อก ๐๓๑๐(๑)/ ๑๖๑๖๘

ลงวันที่ ๒๐ พฤศจิกายน ๒๕๖๖

ก. ผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๖ ราย

๑) นางสาวสุพาทกร จันทร์เปล่ง

๒) นางสาวจันทน์ ไกรภักดิ์ ณ นคร

๓) นายศุภพัทธ์ จิตราชนันท์

๔) นางสาวกนกกร เอบก

๕) นายสุริยา สอนแก้ว

๖) นายวิชาญ ขุนหวัด

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๐๑

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๐๒

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๐๓

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๐๔

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๐๕

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๐๖

สิ่งที่ส่งมาด้วย ๒

เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด เลขทะเบียน ๖-๒๐๔

ที่ อก ๐๓๑๐(๑)/ ๑๖๑๖๘

ลงวันที่ ๒๐ พฤศจิกายน ๒๕๖๖

ข. เจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๑๘๑ ราย

๑) นายภาณุวัฒน์ กิตติคุณาณีย์

๒) นายภัทรพล สว่างใจธรรม

๓) นายธนวิทย์ เทือกคำ

๔) นายศิริโชค พงษ์ประสม

๕) นายณัฐวุฒิ คังแพง

๖) นางสาวจินดา ไชยธรรม

๗) นางสาวสวดี น้อยเจริญ

๘) นางสาวณัฐภาณุ อิมขม

๙) นางสาวนรินทร์ สายแสง

๑๐) นางสาวนันทิ์ สมบูรณ์

๑๑) นางสาวศรัณยา เกลิมอ้างค์

๑๒) นางสาวณัฐพร มงคลจิตร

๑๓) นางสาวศิริลักษณ์ บุญนาค

๑๔) นายณพพงศ์ จันทน์พันธุ์

๑๕) นายณนเศรษฐ์ โทมะสิทธิ์

๑๖) นายธินา จริยา

๑๗) นางสาวภาณุพร แก้วมื่น

๑๘) นางสาวสุวิมล ชัยเรืองวุฒิ

๑๙) นางสาวสุชาดา ธรรมถาวร

๒๐) นางสาวเมธิกา ชัยเดชอนกุล

๒๑) นางสาวศศิธร หนูสวัสดิ์

๒๒) นางสาวเสาวลักษณ์ ภูมิกำพร

๒๓) นายอภิสิทธิ์ สิงหา

๒๔) นายศักดิ์สิทธิ์ โพธิ์พิสุทธิ์

๒๕) ว่าที่ร้อยตรีหญิง พรรณิภา จำเจริญ

๒๖) นางจิตตา คำมูลแก้ว

๒๗) นางสาวอรุณวรรณ รักษ์

๒๘) นางสาวนพรัตน์ แยมกราม

๒๙) นายจุลเดช วารินทร์

๓๐) นางสาวศุภาวิณี รื่องคำ

๓๑) นายพรมณ์ ศรีปัดเนตร

๓๒) นายชุตติ์ สุนธิ์

๓๓) ว่าที่ร้อยตรี เกลิมเกียรติ อมรรศรีเสริม

๓๔) นางสาววริยา สว่างนา

๓๕) นายอนุพงศ์ รัตนศรีประเสริฐ

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๐๑

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๐๒

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๐๓

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๐๔

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๐๕

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๐๖

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๐๗

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๐๘

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๐๙

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๑๐

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๑๑

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๑๒

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๑๓

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๑๔

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๑๕

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๑๖

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๑๗

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๑๘

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๑๙

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๒๐

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๒๑

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๒๒

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๒๓

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๒๔

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๒๕

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๒๖

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๒๗

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๒๘

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๒๙

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๓๐

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๓๑

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๓๒

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๓๓

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๓๔

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๓๕

๑๖) นางสาวจุฑารัตน์...

- ๒ -

๑๖) นางสาวจุฑารัตน์ โอนสันเพียร

๑๗) นางสาวจุฑารัตน์ พิมพ์กิจดิยา

๑๘) นางสาวปรารถนา ศิริกิจดิยา

๑๙) นางสาวศศิธร หนูสวัสดิ์

๒๐) นางสาวศศิธร หนูสวัสดิ์

๒๑) นางสาวศศิธร หนูสวัสดิ์

๒๒) นางสาวศศิธร หนูสวัสดิ์

๒๓) นางสาวศศิธร หนูสวัสดิ์

๒๔) นางสาวศศิธร หนูสวัสดิ์

๒๕) นางสาวศศิธร หนูสวัสดิ์

๒๖) นางสาวศศิธร หนูสวัสดิ์

๒๗) นางสาวศศิธร หนูสวัสดิ์

๒๘) นางสาวศศิธร หนูสวัสดิ์

๒๙) นางสาวศศิธร หนูสวัสดิ์

๓๐) นางสาวศศิธร หนูสวัสดิ์

๓๑) นางสาวศศิธร หนูสวัสดิ์

๓๒) นางสาวศศิธร หนูสวัสดิ์

๓๓) นางสาวศศิธร หนูสวัสดิ์

๓๔) นางสาวศศิธร หนูสวัสดิ์

๓๕) นางสาวศศิธร หนูสวัสดิ์

๓๖) นางสาวศศิธร หนูสวัสดิ์

๓๗) นางสาวศศิธร หนูสวัสดิ์

๓๘) นางสาวศศิธร หนูสวัสดิ์

๓๙) นางสาวศศิธร หนูสวัสดิ์

๔๐) นางสาวศศิธร หนูสวัสดิ์

๔๑) นางสาวศศิธร หนูสวัสดิ์

๔๒) นางสาวศศิธร หนูสวัสดิ์

๔๓) นางสาวศศิธร หนูสวัสดิ์

๔๔) นางสาวศศิธร หนูสวัสดิ์

๔๕) นางสาวศศิธร หนูสวัสดิ์

๔๖) นางสาวศศิธร หนูสวัสดิ์

๔๗) นางสาวศศิธร หนูสวัสดิ์

๔๘) นางสาวศศิธร หนูสวัสดิ์

๔๙) นางสาวศศิธร หนูสวัสดิ์

๕๐) นางสาวศศิธร หนูสวัสดิ์

๕๑) นางสาวศศิธร หนูสวัสดิ์

๕๒) นางสาวศศิธร หนูสวัสดิ์

๕๓) นางสาวศศิธร หนูสวัสดิ์

๕๔) นางสาวศศิธร หนูสวัสดิ์

๕๕) นางสาวศศิธร หนูสวัสดิ์

๕๖) นางสาวศศิธร หนูสวัสดิ์

๕๗) นางสาวศศิธร หนูสวัสดิ์

๕๘) นางสาวศศิธร หนูสวัสดิ์

๕๙) นางสาวศศิธร หนูสวัสดิ์

๖๐) นางสาวศศิธร หนูสวัสดิ์

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๑๖

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๑๗

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๑๘

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๑๙

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๒๐

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๒๑

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๒๒

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๒๓

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๒๔

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๒๕

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๒๖

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๒๗

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๒๘

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๒๙

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๓๐

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๓๑

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๓๒

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๓๓

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๓๔

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๓๕

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๓๖

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๓๗

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๓๘

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๓๙

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๔๐

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๔๑

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๔๒

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๔๓

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๔๔

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๔๕

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๔๖

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๔๗

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๔๘

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๔๙

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๕๐

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๕๑

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๕๒

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๕๓

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๕๔

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๕๕

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๕๖

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๕๗

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๕๘

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๕๙

ทะเบียนเลขที่ ๖-๒๐๔-๓-๐๐๖๐

๑๕) นายประเสริฐ...







ลำดับที่	สารเคมี	วิธีวิเคราะห์
19	Copper	1) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(4)</sup>
20	Cyanide	Distillation, Colorimetric Method <sup>(4)</sup>
21	2,4'-DDD	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
22	4,4'-DDD	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
23	2,4'-DDE	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
24	4,4'-DDE	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
25	2,4'-DDT	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
26	4,4'-DDT	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
27	Dieldrin	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
28	Endosulfan Sulfate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
29	Endosulfan I	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
30	Endosulfan II	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
31	Endrin	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
32	Endrin Aldehyde	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
33	Formaldehyde	Distillation, Colorimetric Method <sup>(4)</sup>
34	Free Chlorine	1) DPD Ferrous Titrimetric Method <sup>(4)</sup> 2) DPD Colorimetric Method <sup>(4)</sup>
35	Heptachlor	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
36	Heptachlor Epoxide	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
37	Hexavalent Chromium	Colorimetric Method <sup>(4)</sup>
38	3-Hydroxycarbofuran	High-Performance Liquid Chromatographic Method <sup>(4)</sup>
39	Lead	1) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(4)</sup>

40 Manganese...

ลำดับที่	สารเคมี	วิธีวิเคราะห์
40	Manganese	1) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(4)</sup>
41	Mercury	1) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(4)</sup>
42	Methiocarb	High-Performance Liquid Chromatographic Method <sup>(4)</sup>
43	Methoxychlor	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
44	Methomyl	High-Performance Liquid Chromatographic Method <sup>(4)</sup>
45	Nickel	1) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(4)</sup>
46	Oil & Grease	1) Liquid-Liquid, Partition-Gravimetric Method <sup>(4)</sup> 2) Soxhlet Extraction Method <sup>(4)</sup>
47	Oxamyl	High-Performance Liquid Chromatographic Method <sup>(4)</sup>
48	Propoxur	High-Performance Liquid Chromatographic Method <sup>(4)</sup>
49	pH	Electrometric Method <sup>(4)</sup>
50	Phenols	1) Distillation, Chloroform Extraction Method <sup>(4)</sup> 2) Distillation, Direct Photometric Method <sup>(4)</sup>
51	Selenium	1) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(4)</sup>
52	Sulfide	Iodometric Method <sup>(4)</sup>
53	Temperature	Laboratory and Field Methods <sup>(4)</sup>
54	Total Dissolved Solids	Dried at 180 °C <sup>(4)</sup>
55	Total Kjeldahl Nitrogen	Semi-Micro Kjeldahl Method <sup>(4)</sup>
56	Total Phosphorous	Digestion, Colorimetric Method <sup>(4)</sup>
57	Total Suspended Solids	Dried from 103-105 °C <sup>(4)</sup>
58	Toxaphene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
59	Trivalent Chromium	1) Digestion, Inductively Coupled Plasma Method; Colorimetric Method; Calculation <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method; Colorimetric Method; Calculation <sup>(4)</sup>
60	Zinc	1) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(4)</sup>

น้ำโคลน...

น้ำโคลน จำนวน 126 รายการ

ลำดับที่	สารเคมี	วิธีวิเคราะห์
1	Acenaphthene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
2	Acetone	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
3	Aldrin	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
4	Anthracene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
5	Antimony	1) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(4)</sup>
6	Arsenic	1) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(4)</sup>
7	Atrazine	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
8	Barium	1) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(4)</sup>
9	Benz(a)anthracene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
10	Benzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
11	Benzo(b)fluoranthene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
12	Benzo(k)fluoranthene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
13	Benzoic Acid	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
14	Benzo(a)pyrene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
15	Benzo(g,h,i)perylene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
16	Beryllium	1) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(4)</sup>
17	Bis(2-chloroethyl)ether	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>

18 Bis(2-ethylhexyl)phthalate...

ลำดับที่	สารเคมี	วิธีวิเคราะห์
18	Bis(2-ethylhexyl)phthalate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
19	Bromodichloromethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
20	Bromoform	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
21	Butanol	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
22	Butyl benzyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
23	Cadmium	1) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(4)</sup>
24	Carbazole	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
25	Carbon disulfide	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
26	Carbon tetrachloride	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
27	Chlordane	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
28	p-Chloroaniline	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
29	Chlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
30	Chlorodibromomethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
31	Chloroform	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
32	2-Chlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
33	Chromium	1) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(4)</sup>
34	Chromium (III)	1) Digestion, Inductively Coupled Plasma Method; Colorimetric Method; Calculation <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method; Colorimetric Method; Calculation <sup>(4)</sup>
35	Chromium (VI)	Colorimetric Method <sup>(4)</sup>

36 Chrysene...



ลำดับที่	สารเคมี	วิธีวิเคราะห์
36	Chrysene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
37	Cyanide	Distillation, Colorimetric Method <sup>(4)</sup>
38	2,4-D	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
39	DDD	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
40	DDE	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
41	DDT	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
42	Dibenz(a,h)anthracene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
43	Di-n-Butyl Phthalate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
44	1,2-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
45	1,3-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
46	1,4-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
47	3,3-Dichlorobenzidine	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
48	1,1-Dichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
49	1,2-Dichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
50	1,1-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
51	cis-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
52	trans-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
53	2,4-Dichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
54	1,2-Dichloropropane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
55	1,3-Dichloropropane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>

56 1,3-Dichloropropene...

ลำดับที่	สารเคมี	วิธีวิเคราะห์
56	1,3-Dichloropropene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
57	Dieldrin	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
58	Diethyl Phthalate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
59	2,4-Dimethylphenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
60	2,4-Dinitrophenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
61	2,4-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
62	2,6-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
63	Di-n-octyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
64	Endosulfan	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
65	Endrin	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
66	Ethylbenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
67	Fluoranthene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
68	Fluorene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
69	Heptachlor	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
70	Heptachlor epoxide	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
71	Hexachlorobenzene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
72	Hexachloro-1,3-butadiene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
73	n-Hexane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
74	α-HCH	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
75	β-HCH	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>

76 γ-HCH...

ลำดับที่	สารเคมี	วิธีวิเคราะห์
76	γ-HCH	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
77	Hexachlorocyclopentadiene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
78	Hexachloroethane	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
79	Indeno(1,2,3-cd)pyrene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
80	Isophorone	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
81	Lead	1) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(4)</sup>
82	Manganese	1) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(4)</sup>
83	Mercury	1) Digestion, Cold Vapor Atomic Absorption Spectrometric Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(4)</sup>
84	Methanol	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
85	Methoxychlor	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
86	Methyl bromide	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
87	Methylene chloride	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
88	2-Methylphenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
89	2-Methylnaphthalene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
90	Methyl tert-butyl Ether	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
91	Naphthalene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
92	Nickel	1) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(4)</sup>
93	Nitrobenzene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>

94 N-Nitrosodiphenylamine...

ลำดับที่	สารเคมี	วิธีวิเคราะห์
94	N-Nitrosodiphenylamine	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
95	N-Nitrosodi-n-Propylamine	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
96	Polychlorinated Biphenyls - PCB 1016 - PCB 1221 - PCB 1232 - PCB 1242 - PCB 1248 - PCB 1254 - PCB 1254	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
97	Pentachlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
98	pH	Electrometric Method <sup>(4)</sup>
99	Phenanthrene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
100	Phenol	1) Distillation, Chloroform Extraction Method <sup>(4)</sup> 2) Distillation, Direct Photometric Method <sup>(4)</sup> 3) Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
101	Pyrene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
102	Selenium	1) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(4)</sup>
103	Silver	1) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(4)</sup>
104	Styrene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
105	1,1,2,2-Tetrachloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
106	Tetrachloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
107	Toluene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
108	Toxaphene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(4)</sup>
109	TPH (C <sub>9</sub> -C <sub>16</sub> )	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(4),23)</sup>

110 TPH (C<sub>9</sub>-C<sub>16</sub>)...



ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
110	TPH (C <sub>8</sub> -C <sub>16</sub> )	Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(9,22)</sup>
111	TPH (C <sub>10</sub> -C <sub>25</sub> )	Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method <sup>(9,22)</sup>
112	1,2,4-Trichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
113	1,1,1-Trichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
114	1,1,2-Trichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
115	Trichloroethylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
116	2,4,5-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
117	2,4,6-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
118	1,3,5-Trimethylbenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
119	Vanadium	1) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(4)</sup>
120	Vinyl acetate	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
121	Vinyl chloride	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
122	m-Xylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
123	o-Xylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
124	p-Xylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
125	Xylene (Total)	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(4)</sup>
126	Zinc	1) Digestion, Inductively Coupled Plasma Method <sup>(4)</sup> 2) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(4)</sup>

อากาศเสีย...

## อากาศเสีย (ปล่องระบาย) จำนวน 28 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Antimony	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(3)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(3)</sup>
2	Arsenic	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(3)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(3)</sup>
3	Beryllium	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(3)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(3)</sup>
4	Cadmium	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(3)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(3)</sup>
5	Carbon Monoxide	1) Instrumental Analyzer Method <sup>(3)</sup> 2) Sampling Bag Non-Dispersive Infrared Method <sup>(3)</sup>
6	Chlorine	1) Absorption Sampling, Ion Chromatographic Method <sup>(3)</sup> 2) Isokinetic Sampling, Ion Chromatographic Method <sup>(3)</sup>
7	Chromium	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(3)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(3)</sup>
8	Cobalt	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(3)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(3)</sup>
9	Copper	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(3)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(3)</sup>
10	Cresol	Absorption Sampling, Gas Chromatographic Method <sup>(3)</sup>
11	Dioxins	Isokinetic Sampling <sup>(3)</sup>
12	Hydrogen Chloride	1) Absorption Sampling, Ion Chromatographic Method <sup>(3)</sup> 2) Isokinetic Sampling, Ion Chromatographic Method <sup>(3)</sup>
13	Hydrogen Fluoride	1) Absorption Sampling, Ion Chromatographic Method <sup>(3)</sup> 2) Isokinetic Sampling, Ion Chromatographic Method <sup>(3)</sup>
14	Hydrogen Sulfide	Absorption Sampling, Iodometric Method <sup>(3)</sup>

15 Lead...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
15	Lead	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(3)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(3)</sup>
16	Manganese	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(3)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(3)</sup>
17	Mercury	1) Isokinetic Sampling, Digestion, Cold-Vapor Atomic Absorption Spectrometric Method <sup>(3)</sup> 2) Isokinetic Sampling, Digestion, Cold-Vapor Atomic Fluorescence Spectrometric Method <sup>(3)</sup>
18	Nickel	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(3)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(3)</sup>
19	Opacity	Ringelmann's Method <sup>(2)</sup>
20	Oxides of Nitrogen	1) Absorption Sampling, Phenoldisulfonic Acid Method <sup>(3)</sup> 2) Absorption Sampling, Alkaline Permanganate/Colorimetric Method <sup>(3)</sup> 3) Instrumental Analyzer Method <sup>(3)</sup>
21	Selenium	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(3)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(3)</sup>
22	Sulfur Dioxide	1) Absorption Sampling, Barium-Thorin Titrimetric Method <sup>(3)</sup> 2) Instrumental Analyzer Method <sup>(3)</sup>
23	Sulfuric Acid	Isokinetic Sampling, Barium-Thorin Titrimetric Method <sup>(3)</sup>
24	Tellurium	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(3)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(3)</sup>
25	Tin	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(3)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(3)</sup>
26	Total Suspended Particulate	1) Isokinetic Sampling, Gravimetric Method <sup>(3)</sup> 2) Paired Train, Isokinetic Sampling, Gravimetric Method <sup>(3)</sup>

27 Vanadium...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
27	Vanadium	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method <sup>(3)</sup> 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(3)</sup>
28	Xylene	Absorption Sampling, Gas Chromatographic Method <sup>(3)</sup>

## สิ่งปฏิกูลหรือวัสดุที่ไม่ใช้แล้ว จำนวน 35 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Aldrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(1,9,24)</sup> 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(9,24)</sup> 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(1,24)</sup>
2	Antimony	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(1,4,16)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(1,4,17)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(7,17)</sup>
3	Arsenic	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(1,4,16)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(1,4,17)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(7,17)</sup>
4	Barium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(1,4,16)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(1,4,17)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(7,17)</sup>

5 Beryllium...



ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
10	Chromium (VI)	1) Waste Extraction, Colorimetric Method <sup>(7,6,15)</sup> 2) Alkaline Digestion, Colorimetric Method <sup>(6,15)</sup>
11	Cobalt	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(1,6,16)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(1,6,17)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(7,17)</sup>
12	Copper	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(1,6,16)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(1,6,17)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(7,17)</sup>
13	2,4-D	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(1,9,26)</sup> 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup> 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(11,26)</sup>
14	DDD	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(1,9,26)</sup> 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup> 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(11,26)</sup>
15	DDE	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(1,9,26)</sup> 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup> 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(11,26)</sup>
16	DDT	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(1,9,26)</sup>

2) Soxhlet..

ลำดับที่	สารพิษ	วิธีการหา
22	Mercury	1) Waste Extraction, Digestion, Cold-Vapor Atomic Absorption Spectrometric Method <sup>(1,6,20)</sup> 2) Waste Extraction, Digestion, Cold-Vapor Atomic Fluorescence Spectrometric Method <sup>(1,6,20)</sup> 3) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method <sup>(20)</sup> 4) Digestion, Cold-Vapor Atomic Fluorescence Spectrometric Method <sup>(20)</sup> 5) Thermal Decomposition Amalgamation and Atomic Absorption Spectrometric Method <sup>(21)</sup>
23	Methoxychlor	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(3,9,20)</sup> 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,20)</sup> 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(1,1,20)</sup>
24	Mirex	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(3,9,20)</sup> 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,20)</sup> 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(1,20)</sup>
25	Molybdenum	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(1,6,16)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(1,6,17)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(7,17)</sup>
26	Nickel	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(3,6,16)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(1,6,17)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(7,17)</sup>
27	Polychlorinated biphenyls (PCBs) - Aroclor 1016 - Aroclor 1221 - Aroclor 1232 - Aroclor 1242 - Aroclor 1248 - Aroclor 1254 - Aroclor 1260	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(1,9,20)</sup> 2) Soxhlet Extraction, Gas Chromatographic Method <sup>(10,20)</sup> 3) Automated Soxhlet Extraction, Gas Chromatographic Method <sup>(3,10)</sup>

- 2-ChlorobiphenylL...



ลำดับที่	สารเคมี	วิธีวิเคราะห์
28	- 2-Chlorobiphenyl - 2,3-Dichlorobiphenyl - 2,2',5-Trichlorobiphenyl - 2,4',5-Trichlorobiphenyl - 2,2',3,5'-Tetrachlorobiphenyl - 2,2',5,5'-Tetrachlorobiphenyl - 2,3',4,4'-Tetrachlorobiphenyl - 2,2',3,4,5'-Pentachlorobiphenyl - 2,2',4,5,5'-Pentachlorobiphenyl - 2,3',4',6-Pentachlorobiphenyl - 2,2',3,4,4',5'-Hexachlorobiphenyl - 2,2',3,4,5,5'-Hexachlorobiphenyl - 2,2',3,5,5',6-Hexachlorobiphenyl - 2,2',4,4',5,5'-Hexachlorobiphenyl - 2,2',3,3',4,4',5-Heptachlorobiphenyl - 2,2',3,4,4',5,5'-Heptachlorobiphenyl - 2,2',3,4,4',5',6-Heptachlorobiphenyl - 2,2',3,4',5,5',6-Heptachlorobiphenyl - 2,2',3,3',4,4',5,5',6-Nonachlorobiphenyl Pentachlorophenol	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(1,9,24)</sup> 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup> 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(11,26)</sup> Electrometric Method <sup>(23,24)</sup> 1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(1,6,16)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(1,6,17)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(7,17)</sup>
29	pH	
30	Selenium	

31 Silver...

ลำดับที่	สารเคมี	วิธีวิเคราะห์
31	Silver	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(1,6,16)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(1,6,17)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(7,17)</sup>
32	Thallium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(1,6,16)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(1,6,17)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(7,17)</sup>
33	Toxaphene	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(9,26)</sup> 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup> 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(11,26)</sup>
34	Vanadium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(1,6,16)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(1,6,17)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(7,17)</sup>
35	Zinc	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method <sup>(1,6,16)</sup> 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(1,6,17)</sup> 3) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(7,17)</sup>

31...

## สินค้าจำนวน 125 รายการ

ลำดับที่	สารเคมี	วิธีวิเคราะห์
1	Acenaphthene	1) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(11,26)</sup>
2	Acetone	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(15,25)</sup> 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method <sup>(13)</sup>
3	Aldrin	1) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(11,26)</sup>
4	Anthracene	1) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(11,26)</sup>
5	Antimony	1) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 2) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(7,17)</sup>
6	Arsenic	1) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 2) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(7,17)</sup>
7	Atrazine	1) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(11,26)</sup>
8	Barium	1) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 2) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(7,17)</sup>
9	Benz(a)anthracene	1) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(11,26)</sup>
10	Benzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(15,25)</sup>

11 Benzo(b)fluoranthene

ลำดับที่	สารเคมี	วิธีวิเคราะห์
11	Benzo(b)fluoranthene	1) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(11,26)</sup>
12	Benzo(k)fluoranthene	1) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(11,26)</sup>
13	Benzoic acid	1) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(11,26)</sup>
14	Benzo(a)pyrene	1) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(11,26)</sup>
15	Benzo(g,h,i)perylene	1) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(11,26)</sup>
16	Beryllium	1) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 2) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method <sup>(7,17)</sup>
17	Bis(2-chloroethyl)ether	1) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(11,26)</sup>
18	Bis(2-ethylhexyl)phthalate	1) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(11,26)</sup>
19	Bromodichloromethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(15,25)</sup>
20	Bromoform	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method <sup>(15,25)</sup>
21	Butanol	Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method <sup>(13,25)</sup>
22	Butyl Benzyl Phthalate	1) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method <sup>(11,26)</sup>

23 Cadmium...



ลำดับที่	สารเคมี	วิธีวิเคราะห์
23	Cadmium	1) Digestion, Inductively Coupled Plasma Method <sup>(7,14)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(7,17)</sup>
24	Carbazole	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
25	Carbon Disulfide	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
26	Carbon tetrachloride	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
27	Chlordane	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
28	p-Chloroaniline	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
29	Chlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
30	Chlorodibromomethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
31	Chloroform	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
32	2-Chlorophenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
33	Chromium	1) Digestion, Inductively Coupled Plasma Method <sup>(7,14)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(7,17)</sup>
34	Chromium (III)	1) Digestion, Inductively Coupled Plasma Method/ Alkaline Digestion, Colorimetric Method; Calculation Method <sup>(7,8,16,19)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method; Alkaline Digestion, Colorimetric Method; Calculation Method <sup>(7,8,17,19)</sup>
35	Chromium (VI)	Alkaline Digestion, Colorimetric Method <sup>(8,19)</sup>

35 Chrysene...

ลำดับที่	สารเคมี	วิธีวิเคราะห์
36	Chrysene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
37	Cyanide	Extraction, Distillation, Colorimetric Method <sup>(27,28,29)</sup>
38	2,4-D	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
39	DDD	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
40	DDE	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
41	DDT	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
42	Dibenz(a,h)anthracene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
43	Di-n-Butyl Phthalate	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
44	1,2-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
45	1,3-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
46	1,4-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
47	3,3-Dichlorobenzidine	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
48	1,1-Dichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>

49 1,2-Dichloroethane...

ลำดับที่	สารเคมี	วิธีวิเคราะห์
49	1,2-Dichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
50	1,1-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
51	cis-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
52	trans-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
53	2,4-Dichlorophenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
54	1,2-Dichloropropane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
55	1,3-Dichloropropane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
56	1,3-Dichloropropene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
57	Dieldrin	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
58	Diethyl Phthalate	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
59	2,4-Dimethylphenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
60	2,4-Dinitrophenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
61	2,4-Dinitrotoluene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
62	2,6-Dinitrotoluene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>

63 Di-n-Octyl Phthalate...

ลำดับที่	สารเคมี	วิธีวิเคราะห์
63	Di-n-Octyl Phthalate	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
64	Endosulfan	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
65	Endrin	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
66	Ethylbenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
67	Fluoranthene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
68	Fluorene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
69	Heptachlor	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
70	Heptachlor epoxide	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
71	Hexachlorobenzene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
72	Hexachloro-1,3-butadiene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
73	n-Hexane	1) Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup> 2) Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method <sup>(19)</sup>

73 n-Hexane...



ลำดับที่	สารเคมี	วิธีวิเคราะห์
74	$\alpha$ -HCH	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
75	$\beta$ -HCH	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
76	$\gamma$ -HCH	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
77	Hexachlorocyclopentadiene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
78	Hexachloroethane	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
79	Indeno(1,2,3-cd)pyrene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
80	Isophorone	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
81	Lead	1) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(7,17)</sup>
82	Manganese	1) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(7,17)</sup>
83	Mercury	1) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method <sup>(20)</sup> 2) Thermal Decomposition, Amalgamation, and Atomic Absorption Spectrophotometry <sup>(21)</sup> 3) Digestion, Cold-Vapor Atomic Fluorescence Spectrometric Method <sup>(20)</sup>

84 Methanol...

ลำดับที่	สารเคมี	วิธีวิเคราะห์
84	Methanol	1) Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup> 2) Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method <sup>(13,25)</sup>
85	Methoxychlor	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
86	Methyl Bromide	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
87	Methylene Chloride	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
88	2-methylphenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
89	2-Methylnaphthalene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
90	Methyl, tert-Butyl Ether	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
91	Naphthalene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
92	Nickel	1) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(7,17)</sup>
93	Nitrobenzene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
94	N-Nitrosodiphenylamine	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
95	N-Nitrosodi-n-propylamine	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>

96 Polychlorinated biphenyls (PCBs)

ลำดับที่	สารเคมี	วิธีวิเคราะห์
96	Polychlorinated biphenyls (PCBs) - Aroclor 1016 - Aroclor 1221 - Aroclor 1232 - Aroclor 1242 - Aroclor 1248 - Aroclor 1254 - Aroclor 1260 - 2-Chlorobiphenyl - 2,2',3,5'-Tetrachlorobiphenyl - 2,2',5,5'-Tetrachlorobiphenyl - 2,3',4,4'-Tetrachlorobiphenyl - 2,2',3,4,5'-Pentachlorobiphenyl - 2,2',4,5,5'-Pentachlorobiphenyl - 2,3',4',6-Pentachlorobiphenyl - 2,2',3,4,4',5'-Hexachlorobiphenyl - 2,2',3,4,5,5'-Hexachlorobiphenyl - 2,2',3,5,5',6-Hexachlorobiphenyl - 2,2',4,4',5,5'-Hexachlorobiphenyl - 2,2',3,3',4,4',5-Heptachlorobiphenyl - 2,2',3,4,4',5,5'-Heptachlorobiphenyl - 2,2',3,4,4',5,6-Heptachlorobiphenyl - 2,2',3,4',5,5',6-Heptachlorobiphenyl - 2,2',3,3',4,4',5,6-Nonachlorobiphenyl	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
97	Pentachlorophenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
98	Phenanthrene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>

99 Phenol...

ลำดับที่	สารเคมี	วิธีวิเคราะห์
99	Phenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
100	Pyrene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
101	Selenium	1) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(7,17)</sup>
102	Silver	1) Digestion, Inductively Coupled Plasma Method <sup>(7,16)</sup> 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method <sup>(7,17)</sup>
103	Styrene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
104	1,1,2,2-Tetrachloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
105	Tetrachloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
106	Toluene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
107	Toxaphene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(10,26)</sup> 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method <sup>(11,26)</sup>
108	TPH (C <sub>6</sub> -C <sub>8</sub> )	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
109	TPH (C <sub>9</sub> -C <sub>16</sub> )	1) Automate Extraction, Gas Chromatographic Method <sup>(12,22)</sup> 2) Solvent Extraction, Gas Chromatographic Method <sup>(12,22)</sup> 3) Ultrasonic Extraction, Gas Chromatographic Method <sup>(12,23)</sup>
110	TPH (C <sub>16</sub> -C <sub>35</sub> )	1) Automate Extraction, Gas Chromatographic Method <sup>(12,22)</sup> 2) Solvent Extraction, Gas Chromatographic Method <sup>(12,22)</sup> 3) Ultrasonic Extraction, Gas Chromatographic Method <sup>(12,23)</sup>
111	1,2,4-Trichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
112	1,1,1-Trichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
113	1,1,2-Trichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>
114	Trichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method <sup>(15,25)</sup>

115 2,4,5-Trichlorophenol...





อนึ่ง หนังสือฉบับนี้จะหมดอายุพร้อมหนังสือต่ออายุฉบับใหม่เมื่อพ้นกำหนดปฏิบัติภารกิจพิเศษ  
ในวันที่ ๒ กันยายน ๒๕๖๔

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

  
(นายพรศ ก่อกรอง)  
รองอธิบดี ปฏิบัติราชการแทน  
อธิบดีกรมโรงงานอุตสาหกรรม

กองวิจัยและเตือนภัยมลพิษโรงงาน  
กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ  
โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕  
โทรสาร ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๑๔๙  
ไปรษณีย์อิเล็กทรอนิกส์ saraban@dlw.mail.go.th



ที่ อก ๐๓๓๐(๑)/ ๑๒๓๖ ๘ /

กรมโรงงานอุตสาหกรรม  
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท  
เขตราชเทวี กรุงเทพฯ ๑๐๕๐๐

๑๘ ธันวาคม ๒๕๖๓

เรื่อง ยกเลิกบุคลากรของห้องปฏิบัติการวิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท เอแอลเอส แลบริทอรี่ กรุ๊ป (ประเทศไทย) จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน  
ลงวันที่ ๒ ธันวาคม ๒๕๖๓


ตามคำขอที่อ้างถึง บริษัท เอแอลเอส แลบริทอรี่ กรุ๊ป (ประเทศไทย) จำกัด ห้องปฏิบัติการ  
วิเคราะห์เอกชน เลขทะเบียน ๖-๒๐๔ สถานที่ตั้งเลขที่ ๑๐๔ ซอยพัฒนาการ ๔๐ ถนนพัฒนาการ แขวงพัฒนาการ  
เขตสวนหลวง กรุงเทพมหานคร ขอยกเลิกบุคลากร ความละเอียดแจ้งแล้ว นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้ออกเลิกเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์  
จำนวน ๘ ราย ได้แก่

๑) นายประพนธ์ วรรณชัย	ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๖๐
๒) นายจิรเมธ ขวาละออ	ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๗๒
๓) นายพิรพัฒน์ กำคำ	ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๐๘
๔) นางสาวอรยา คำค่อง	ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๓๔
๕) นายกิตติพงศ์ แซ่ลี	ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๔๔
๖) นายจิรเมธ ประเสริฐศิริพงศ์	ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๖๐
๗) นายภัทรพงษ์ มณฑาทอง	ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๖๗
๘) นางสาวจาวรรณ กระจำพันธุ์	ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๘๑

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

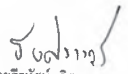  
(นายธีรทัศน์ อิศรางกูร ณ อยุธยา)  
รองอธิบดี ปฏิบัติราชการแทน  
อธิบดีกรมโรงงานอุตสาหกรรม

กองวิจัยและเตือนภัยมลพิษโรงงาน  
กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ  
โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕  
โทรสาร ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๑๔๙  
ไปรษณีย์อิเล็กทรอนิกส์ saraban@dlw.mail.go.th



"อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว"



"อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว"

