

ภาคผนวก ง

ใบรับรองการสอบเทียบเครื่องมือ



right solutions.
right partner.

รายการเครื่องมือที่ใช้ในการวิเคราะห์ / ทดสอบ

Sample Name	Parameter	Equipment Name	ID No.	Calibrated Date	Next Cal	Freq. Calibrate (Months)
Songkhla Lab	Residual Free Chlorine	Chlorine Meter	PHK_FS0003	13-Aug-24	13-Aug-25	12
Songkhla Lab	Iron	ICP-MS	SGK_CL0048	2-Aug-23	2-Feb-25	18
Songkhla Lab	Iron	Cold Room Water	SGK_CL0065	1-Jul-24	1-Jan-26	18
Songkhla Lab	Conductivity	Conductivity Meter	SGK_CL0032	6-May-24	6-May-25	12
Songkhla Lab	pH at 25 °C	pH meter	SGK_CL0030	19-Oct-24	19-Apr-26	18
Songkhla Lab	Total Dissolved Solids 180 °C	Electronic Top-Loading Balance	SGK_CL0045	15-Jan-24	15-Jan-25	12
Songkhla Lab	Total Dissolved Solids 180 °C	Oven	SGK_CL0024	19-Oct-24	19-Apr-26	18
Songkhla Lab	Total Dissolved Solids 103-105 °C	Electronic Top-Loading Balance	SGK_CL0045	15-Jan-24	15-Jan-25	12
Songkhla Lab	Total Dissolved Solids 103-105 °C	Oven	SGK_CL0024	19-Oct-24	19-Apr-26	18
Songkhla Lab	Turbidity	Turbidity Meter	SGK_FS0045	22-May-24	22-May-25	12
Songkhla Lab	BOD	Incubator	SGK_CL0028	13-Jul-23	13-Jan-25	18
Songkhla Lab	BOD	DO/BOD Analyser	SGK_CL0073	21-May-24	21-Nov-25	18
Songkhla Lab	Oil & Grease	Electronic Top-Loading Balance	SGK_CL0045	15-Jan-24	15-Jan-25	12
Songkhla Lab	Oil & Grease	Oven	SGK_CL0024	19-Oct-24	19-Apr-26	18
Songkhla Lab	Oil & Grease	Water Bath	SGK_CL0035	13-Jul-23	13-Jan-25	18
Songkhla Lab	Sulfide	Cold Room Water	SGK_CL0065	1-Jul-24	1-Jan-26	18
Songkhla Lab	Total Suspended Solids	Electronic Top-Loading Balance	SGK_CL0045	15-Jan-24	15-Jan-25	12
Songkhla Lab	Total Suspended Solids	Oven	SGK_CL0024	19-Oct-24	19-Apr-26	18
Songkhla Lab	Mixed Liquor Suspended Solids	Oven	SGK_CL0024	19-Oct-24	19-Apr-26	18
Songkhla Lab	Total Coliform	Autoclave	SGK_ML0001	1-Jul-24	1-Jan-26	18
Songkhla Lab	Total Coliform	Incubator	SGK_ML0013	15-Jan-24	15-Jul-25	18
Songkhla Lab	Total Coliform	pH Meter	SGK_ML0016	2-Jul-24	2-Jan-26	18
Songkhla Lab	Total Coliform	Water Bath	SGK_ML0021	1-Jul-24	1-Jan-26	18
Songkhla Lab	<i>Escherichia coli</i>	Autoclave	SGK_ML0002	1-Jul-24	1-Jan-26	18
Songkhla Lab	<i>Escherichia coli</i>	Incubator	SGK_ML0013	15-Jan-24	15-Jul-25	18
Songkhla Lab	<i>Escherichia coli</i>	pH meter	SGK_ML0016	2-Jul-24	2-Jan-26	18
Songkhla Lab	<i>Escherichia coli</i>	Water Bath	SGK_ML0021	1-Jul-24	1-Jan-26	18
Songkhla Lab	Standard Plate Count	Incubator	SGK_ML0013	15-Jan-24	15-Jul-25	18
Water Lab	Bicarbonate Alkalinity	Burette	BKK_EN0171	27-Feb-24	27-Aug-25	18
Water Lab	Carbonate Alkalinity	Burette	BKK_EN0171	27-Feb-24	27-Aug-25	18
Water Lab	Chloride	Ion Chromatography	BKK_EN0069	12-Jan-24	12-Jan-25	12
Water Lab	Hydroxide Alkalinity	Burette	BKK_EN0171	27-Feb-24	27-Aug-25	18
Water Lab	Phenolphthalein Alkalinity	Burette	BKK_EN0171	27-Feb-24	27-Aug-25	18
Water Lab	Methyl Orange Alkalinity	Burette	BKK_EN0171	27-Feb-24	27-Aug-25	18
Water Lab	Residual Alkalinity	Burette	BKK_EN0171	27-Feb-24	27-Aug-25	18
Water Lab	Permanent Hardness	Burette	BKK_EN0171	27-Feb-24	27-Aug-25	18
Water Lab	Nitrate	Ion Chromatography	BKK_EN0069	12-Jan-24	12-Jan-25	12
Water Lab	Sulfate	Ion Chromatography	BKK_EN0069	12-Jan-24	12-Jan-25	12
Water Lab	Temporary Hardness	Burette	BKK_EN0171	27-Feb-24	27-Aug-25	18
Water Lab	Total Alkalinity	Burette	BKK_EN0171	27-Feb-24	27-Aug-25	18
Water Lab	Total Hardness	Burette	BKK_EN0171	27-Feb-24	27-Aug-25	18
Water Lab	Total Kjeldahl Nitrogen	Digestion Unit	BKK_EN0223	24-May-24	24-May-25	12
Water Lab	Total Kjeldahl Nitrogen	Discrete analyzer	BKK_EN0037	16-Aug-24	16-Aug-25	12
Water Lab	<i>Legionella spp.</i>	Balance	BKK_ML0001	23-Jan-24	23-Jan-25	12
Water Lab	<i>Legionella spp.</i>	Autoclave	BKK_ML0041	4-Oct-23	4-Apr-25	18
Water Lab	<i>Legionella spp.</i>	Incubator	BKK_ML0189	26-Sep-24	26-Mar-26	18
Water Lab	<i>Legionella spp.</i>	Hot Air Oven	BKK_ML0013	23-Apr-24	23-Oct-25	18

HACH COMPANY

C/O AB Sciex (Thailand) Limited, Building D Room No. D3 11, 3rd Floor , No. 735/4, Srinakarin Road, Pattanakarn, Suanluang, Bangkok
 [Phone +66 (02) 026-3529 Ext. 0 | Fax +66(02) 026-3572 | www.sea.hach.com]

LABX 2302163

Test Report

Customers	:	ALS Laboratory Group (Thailand) Co., Ltd.			
Equipment	:	Chlorine Meter	Manufacturer	:	HACH
Controller Model	:	DR300	ID No.	:	PHK_FS0003
Controller Serial No.	:	19050A001118	Sensor Serial No.	:	-
Date of test	:	07/08/2023	Period	:	-
Environment temperature	:	25.0 °C	Humidity	:	56.0 %RH

Results

Instrument Checked


Item	Characteristic	Before	After	Remark
1	Visual Inspect	☒ Pass ☐ Fail	☒ Pass ☐ Fail	
2	Power Supply (4.5 – 6.0 VDC)	6.0 VDC	6.0 VDC	
3	Display Check	☒ Pass ☐ Fail	☒ Pass ☐ Fail	
4	Keyboard Check	☒ Pass ☐ Fail	☒ Pass ☐ Fail	
5	Function System Program	☒ Pass ☐ Fail	☒ Pass ☐ Fail	

Warning and Error Checked

Item	Event	Before	After
6	Error list	☒ None ☐ Appear	☒ None ☐ Appear

Check with Standard

Item	Characteristic	Before	After	Remark
	DPD-CHLORINE-LR			
7	Blank (0.00 mg/l)	0.00 mg/l	0.00 mg/l	
8	Standard Cl2 No. 1 (0.19 ± 0.09 mg/l)	0.16 mg/l	0.18 mg/l	
9	Standard Cl2 No. 2 (0.87 ± 0.10 mg/l)	0.84 mg/l	0.85 mg/l	
10	Standard Cl2 No. 3 (1.55 ± 0.14 mg/l)	1.53 mg/l	1.54 mg/l	
	DPD-CHLORINE-HR			
11	Blank (0.0 mg/l)	0.0 mg/l	0.0 mg/l	
12	Standard Cl2 No. 1 (2.3 ± 0.2 mg/l)	2.2 mg/l	2.2 mg/l	
13	Standard Cl2 No. 2 (3.9 ± 0.3 mg/l)	3.8 mg/l	3.8 mg/l	
14	Standard Cl2 No. 3 (6.7 ± 0.6 mg/l)	6.7 mg/l	6.7 mg/l	

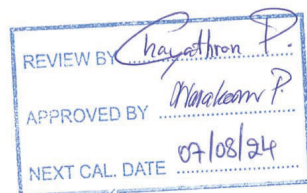

HACH COMPANY

C/O AB Sciex (Thailand) Limited, Building D Room No. D3 11, 3rd Floor , No. 735/4, Srinakarin Road, Pattanakarn, Suanluang, Bangkok
 [Phone +66 (02) 026-3529 Ext. 0 | Fax +66(02) 026-3572 | www.sea.hach.com]

LABX 2302163

Summary of checked

- ☒ The instrument can work normally and efficiently. (เครื่องมือวัดสามารถทำงานได้ปกติและมีประสิทธิภาพ)
- ☐ The instrument can work but it's requiring to maintenance. (เครื่องมือวัดสามารถทำงานได้แต่ต้องบำรุงรักษา)
- ☐ The instrument could not work it's requiring to repair. (เครื่องมือวัดไม่สามารถทำงานได้และต้องการซ่อมบำรุง)

Remark:

Standard Equipment Used

Equipment	Equipment I.D.
Standard Chlorine DPD-CHLORINE-LR	Lot No. A3020 Exp date : Feb-25
Standard Chlorine DPD-CHLORINE-HR	Lot No. A3041 Exp date : Feb-25
Digital multi meter	S/N : 21190066 Due date : Jun-24
Thermo hygrometer	S/N : 45146347 Due date : Aug-24

Test By : WILAILAK S.
 (Miss Wilailak Sawangpun)
 Service Engineer

Approved by : S. Suanun
 (Mr. Suanun Sartyangkool)
 Position : Assistant Service Division Manager





REVIEW BY *Eng. Anthon D*
APPROVED BY *Moraporn P*
NEXT CAL. DATE 18/08/25

Test Report

Customers	ALS Laboratory Group (Thailand) Co., Ltd.		
Equipment	Colorimeter	Manufacturer	HACH
Controller Model	☒ DR300 ☐ Pocket II	ID No.	PHK_FS0003
Controller Serial No.	19050A001118	Sensor Serial No.	-
Date of test	13/08/2024	Period	1 Year
Environment temperature	25 °C	Humidity	60 %RH

Results

Instrument Checked

Item	Characteristic	Before		After		Remark
1	Visual Inspect	☒ Pass	☐ Fail	☒ Pass	☐ Fail	
2	Power Supply (4.5 - 6.0 VDC)	6.0	VDC	6.0	VDC	
3	Display Check	☒ Pass	☐ Fail	☒ Pass	☐ Fail	
4	Keyboard Check	☒ Pass	☐ Fail	☒ Pass	☐ Fail	
5	Function System Program	☒ Pass	☐ Fail	☒ Pass	☐ Fail	

Warning and Error Checked

Item	Event	Before		After	
6	Error list	☒ None	☐ Appear	☒ None	☐ Appear

Check with Standard

Item	Characteristic	Before		After		Remark
DPD-CHLORINE-LR						
7	Blank (0.00 mg/l)	0.00	mg/l	0.00	mg/l	
8	Standard Cl2 No. 1 (0.19 ± 0.09 mg/l)	0.21	mg/l	0.20	mg/l	
9	Standard Cl2 No. 2 (0.87 ± 0.10 mg/l)	0.87	mg/l	0.87	mg/l	
10	Standard Cl2 No. 3 (1.55 ± 0.14 mg/l)	1.55	mg/l	1.55	mg/l	
DPD-CHLORINE-HR						
11	Blank (0.0 mg/l)	0.0	mg/l	0.0	mg/l	
12	Standard Cl2 No. 1 (2.2 ± 0.2 mg/l)	2.2	mg/l	2.2	mg/l	
13	Standard Cl2 No. 2 (3.9 ± 0.3 mg/l)	3.9	mg/l	3.9	mg/l	
14	Standard Cl2 No. 3 (6.8 ± 0.6 mg/l)	6.8	mg/l	6.8	mg/l	



Summary of checked

- ☒ The instrument can work normally and efficiently. (เครื่องมือวัดสามารถทำงานได้ปกติและมีประสิทธิภาพ)
☐ The instrument can work but it's requiring to maintenance. (เครื่องมือวัดสามารถทำงานได้แต่ต้องบำรุงรักษา)
☐ The instrument could not work it's requiring to repair. (เครื่องมือวัดไม่สามารถทำงานได้และต้องการซ่อมบำรุง)

Remark:

Standard Equipment Used

Equipment	Equipment I.D.	
Standard Absorbance DPD-CHLORINE-LR	Lot No. A3020	Exp date : Feb-25
Standard Absorbance DPD-CHLORINE-HR	Lot No. A4025	Exp date : Jan-26
Digital multi meter	S/N : 19430676	Due date : Oct-24
Thermo hygrometer	S/N : 34855953	Due date : Oct-24

Test By :

WILAILAK S.

(Miss Wilailak Sawangpun)

Service Engineer

Approved by :

S. S.

(Mr. Suanun Sartyangkool)

Position : Assistant Service Division Manager



Agilent Technologies (Thailand) Limited
U CHU LIANG BLDG. 22/F UNIT A.D
968 RAMA 4 ROAD, SILOM, BANGRAK
Bangkok 10500 Thailand

Tel. +662 637 6363
Fax: +662 632 4334
Email: ccc-smi@agilent.com
Website: www.agilent.com/chem

Service Confirmation Number: 6905074481

Service Confirmation Date: 02.08.2023

Customer Contact:

ALS Laboratory Group (Thailand) Co
Ltd
Branch Number 0002
114/1 Moo8 Banplu Subdistrict Hat
Yai District
TAX ID : 0105540004859
kanitta.hemprasatpor@alsglobal.com
0811721334

Invoice To:

ALS Laboratory Group (Thailand) Co
Ltd
Branch Number 0002
114/1 Moo8 Banplu Subdistrict Hat
Yai District SONGKHLA 90250

Payer:

ALS Laboratory Group (Thailand) Co
Ltd Head Office
104 Phatthanakan 40 Phatthanakan Rd
Khwaeng Phatthanakan Khet Suan
Luang
BANGKOK 10250

Delivery Site:

ALS Laboratory Group (Thailand) Co
Ltd
Branch Number 0002
114/1 Moo8 Banplu Subdistrict Hat
Yai District

Location:

Room
Bldg
Lab
Dept

SERVICE REPORT

Customer Purchase Order Number:	Customer Number: 70579367
Service Request:	Service Request Date:
Service Order: 6005939422	Service Confirmation: 6905074481

REVIEW BY Kavinda J.
APPROVED BY Kanitta H.
NEXT CAL. DATE 2/02/25

Direct Inquiries to:

Contact Name: Customer Contact Center
Contact E-mail: ccc-smi@agilent.com
Contact Telephone: +662 637 6363
Contact Fax: +662 632 4334

Service Instrument:

Model Number	Model Description	Serial Number	System Handle	Parent Asset
SYS-IM- 7900	ICPMS 7900 System			
G7201C	ICP-MS MassHunter SW only (excludes PC)	USH3789575	ICP MS 7900	SYS-IM-7900
G8403A	Agilent 7900 ICP-MS	JP16511669	ICP MS 7900	SYS-IM-7900
G8411A	ISIS 3 for Agilent 7850/7900/8900	JP16510379	ICP MS 7900	SYS-IM-7900

Service Items:

Item	Service/Part #	Description	Qty	Entitlement	Service Start	Service End
1000	EQO	Enterprise Operational Qualification	1.00	Agreement Entitlement - 100 % covered	02.08.2023	02.08.2023
1010	5185-5850	ICP-MS Checkout Solutions	1.00	Agreement Entitlement - 100 % covered		

Additional Information:

products | applications | software | services

Learn more about Agilent's Special Offers, Products, Services and our
full range of laboratory productivity solutions optimized for your
applications and workflows. Visit us at www.agilent.com/chem

Agilent Technologies (Thailand) Limited, Head Office
U Chu Liang Bldg. 22/F Unit A.D
968 Rama 4 Road, Silom, Bangkok
Bangkok 10500 Thailand
Tax ID : 0105542068218


Citibank N.A. Bangkok Branch
399 Interchange 21 Building, Sukhumvit Road, Klongtoey Nau
Sub-district, Wattana District, Bangkok 10110 Thailand
Acc. No: 012-4452-007,
THB:Krung Thai Bank PCL
Siam Square Br., 416/1-2 Rama I Rd., Pathumwan, BKK 10330
Thailand

ORIGINAL

Service Confirmation Number: 6905074481

Service Confirmation Date: 02.08.2023

Service Information:

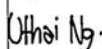Problem Description: WU-S-QQ-IM-7900-5001093854		
Service Provided: -Perform OQ hardware. -Test OQ of instrument ICP-MS = SGK_CL0048. All tests Passed.		
Service Overview Code: Reason Code: Scheduled Service Diagnosis Code: Scheduled Service Resolution Code: Scheduled Service		
Reported Hours: 9.0	Travel Hours: 3.0	
Customer Field Service Representative Name: Uthai Ngamlersirichai	Customer Field Service Representative Signature: 	Date: 02 Aug 2023
Customer Name: KANITTA HEMPRASATPORN	Customer Signature: 	Date: 02 Aug 2023
Additional Comments:		

Southern Calibration Service Co., Ltd.

669/35 Karnjanavanit Rd., Banpru, Hatyai, Songkla 90250 Thailand
Tel : 08 1599 0417 Fax : 0 7480 5133 Email : s.calibration@gmail.com www.scal-lab.com



CALIBRATION CERTIFICATE

Issued Date : 4-Jul-2024

Certificate No. : 24TH2757

CSR No. : A150/07473

Page. : 1 of 3

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Karnchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Chemical Laboratory
Instrument Name : Cold Room
Manufacturer : Danfoss
Model : N/A
Serial No. : N/A
ID No. : SGK_CL0065
Resolution : 0.1 °C
Received Date : 1-Jul-2024
Calibrated Date : 1-Jul-2024
Ambient Temperature : (30 ± 10) °C
Relative Humidity : (50 ± 30) %

REVIEW BY Ananta B.
APPROVED BY Kanitta H.
NEXT CAL. DATE 1/09/26

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WI.012 based on GLA - 20

The Southern Calibration Service Co.,Ltd.calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and /or national standards which realize the units of measurement according to the International System of Unit (SI) through :

- TISTR : Thailand Institute of Scientific and Technological Research

Calibrated by : Ibrahim Saleemin

Approved by :


Imron Rattanayutim / Technical Manager



The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 24TH2757
CSR No. : A150/07473
Page. : 2 of 3

Details of Calibration


1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Data Acquisition/Switch Unit	34970A	MY58009813	PSL-T0707-1/67	22-May-2025

2. The results reported in this certificate refer to the condition of the instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.
3. This certificate is not certified any commercial transaction
4. Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration : (✓) Without Adjustment () After Adjustment

1. Sensor Installation Diagram

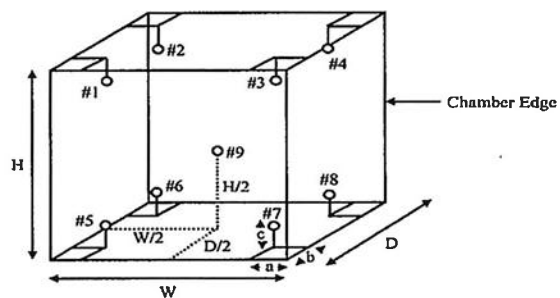

Sensor Installation Details

a = 5.0 cm
b = 5.0 cm
c = 5.0 cm

Dimension of the chamber

W = 40.0 cm
H = 40.0 cm
D = 33.0 cm



Certificate No. : 24TH2757
CSR No. : A150/07473
Page. : 3 of 3

Result of Calibration :

2. Temperature Measurement Accuracy Test

The measurement results of the Cold Room and associates are reported in the manner as shown below

Cal point (°C)	Measured Standard Temperature At Spread Locations (°C)									Uncertainty (± °C)
	#1	#2	#3	#4	#5	#6	#7	#8	Ref. 9	
4	3.90	4.04	3.72	3.78	3.69	3.76	3.55	3.43	3.34	0.39

3. Performance Result

The performance of the Cold Room are reported as shown below

Cal point (°C)	UUC Setting (°C)	UUC Reading (°C)	Temperature Stability (± °C)	Temperature Uniformity (°C)	Overall Variation (°C)
4	4.0	4.0	0.10	0.75	0.75

- UUC = Unit Under Calibration

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%

... End ...



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Cert.No.: 24CH541

Page.: 2 of 2

Certificate of Calibration

Cert.No.: 24CH541

Page.: 1 of 2

Equipment : Conductivity Meter
Manufacturer : Mettler Toledo
Model : SevenCompact S230
Serial No. : B608134488
ID No. : SGK_CL0032
Condition As-Received: Used Item
Received Date : 29 April 2024
Calibration Date : 06 May 2024
Reference : 2404-0708DSC-2
Submitted by : ALS Laboratory Group (Thailand) Co.,Ltd. Songkhla Branch.
114/1 Moo 8 Kamchanawanich Rd.,
T.Ban Phru, A.Hat Yai, Songkhla 90250 Thailand
Ambient Temperature : (25 ± 2.5) °C
Relative Humidity : (50 ± 15) %
Calibration Procedure: In-house method :
- CP-CH6 by direct measurement
with certified reference material (CRM)
Calibrated by : Warakom Lerngagtrakul
Approved by : 
Approved Signatory
() Unnopphol Harachai
() Ponpan Paipim
(✓) Saithip Meangmai
Issue Date : 08 May 2024

REVIEW BY *Kanitta H*
APPROVED BY *Ananta B*
NEXT CAL DATE **06/05/25**

Condition of this result of calibration

1. Reference Standard Instrument :-

Instrument	Serial No.	ID No.	Certificate No.	Due date
1) Thermometer	1963878	130RC095	231051	05 Sep 2024

- This Certification is traceable to SI Through Technology Promotion Association (Thailand - Japan)

2. Certified Reference Materials :-

- Conductivity calibration solution, CPA chem Ltd., The measurement results are traceable to SI through CPA chem Ltd., ANSI-ASQ National Accreditation Board, Accredited No. AR-1835

Conductivity Solution	Manufacturer	Lot No.	Exp. date
84.00 µS/cm	CPA Chem	931954	01 Oct 2024
1413.0 µS/cm	CPA Chem	940111	02 Nov 2024
12.880 mS/cm	CPA Chem	936625	19 Oct 2024

- Control Conductivity calibration solution temperature by Water bath (25 ± 0.1) °C

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration results

Function : Conductivity Measurement

(*) Without Adjustment

Conductivity Electrode Serial No.: 5820120321

Standard Conductivity Solution	UUC* Reading	Uncertainty of Measurement (±)	Coverage factor k
84.00 µS/cm	86.4 µS/cm	0.62 µS/cm	2.00
1413.0 µS/cm	1433 µS/cm	9.2 µS/cm	2.00
12.880 mS/cm	12.88 mS/cm	0.086 mS/cm	2.00

Remark : - UUC* = Unit Under Calibration

- Cell constant = 0.574487 cm⁻¹

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-o0o-

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert.No.: 24CH542

Page.: 1 of 2

Equipment : Conductivity Meter
Manufacturer : Mettler Toledo
Model : SevenCompact S230
Serial No. : B608134488
ID No. : SGK_CL0032
Condition As-Received: Used Item
Received Date : 29 April 2024
Calibration Date : 08 May 2024
Reference : 2404-0708DSC-3
Submitted by : ALS Laboratory Group (Thailand) Co.,Ltd. Songkhla Branch.
114/1 Moo 8 Kamchanawanich Rd.,
T.Ban Phru, A.Hat Yai, Songkhla 90250 Thailand
Ambient Temperature : (25 ± 2.5) °C
Relative Humidity : (50 ± 15) %
Calibration Procedure: In -house method :
- CP-CH6 by direct measurement
with certified reference material (CRM)
Calibrated by : Warakorn Lemgagtrakul
Approved by : 
Approved Signatory
() Unnopphol Harachai
() Ponpan Paipim
(✓) Saithip Meangmai
Issue Date : 08 May 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Cert.No.: 24CH542

Page.: 2 of 2

Condition of this result of calibration

1. Reference Standard Instrument :-

Instrument	Serial No.	ID No.	Certificate No.	Due date
1) Thermometer	1963878	130RC095	231051	05 Sep 2024

- This Certification is traceable to SI Through Technology Promotion Association (Thailand - Japan)

2. Certified Reference Materials :-

- Conductivity calibration solution, CPA chem Ltd., The measurement results are traceable to SI through CPA chem Ltd., ANSI-ASQ National Accreditation Board, Accredited No. AR-1835

Conductivity Solution	Manufacturer	Lot No.	Exp. date
84.00 µS/cm	CPA Chem	931954	01 Oct 2024

- Control Conductivity calibration solution temperature by Water bath (25 ± 0.1) °C

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration results

Function : Conductivity Measurement

(*) After Adjustment at 84.00 µS/cm

Conductivity Electrode Serial No.: 5816340156

Standard Conductivity Solution	Before Adjustment UUC* Reading	After Adjustment UUC* Reading	Uncertainty of Measurement (±)	Coverage factor k
84.00 µS/cm	85.8 µS/cm	84.0 µS/cm	0.62 µS/cm	2.00

Remark : - UUC* = Unit Under Calibration

- Cell constant = 0.097780 cm⁻¹

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-



Southern Calibration Service Co., Ltd.

669/35 Karnjanavanit Rd., Banpru, Hatyai, Songkla 90250 Thailand
Tel : 08 1599 0417 Fax : 0 7480 5133 Email : s.calibration@gmail.com www.scal-lab.com



CALIBRATION CERTIFICATE

Issued Date : 1-May-2023

Certificate No. : 23CH0203

CSR No. : A088/04367

Page : 1 of 2

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Kamchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Chemical Laboratory
Instrument Name : pH meter
Manufacturer : Mettler Toledo
Model : S220
Serial No. : B625631849
ID No. : SGK_CL0030
Electrode No. : 1204613
Received Date : 28-Apr-2023
Calibrated Date : 28-Apr-2023
Ambient Temperature : (25 ± 3) °C
Relative Humidity : (55 ± 15) %

REVIEW BY	Ananta B.
APPROVED BY	Kampha H.
NEXT CAL. DATE	28 / 10 / 2024

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WI.008 based on direct measurement by using certified reference Material (CRM)

The Southern Calibration Service Co., Ltd. calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and /or national standards which realize the units of measurement according to the International System of Unit (SI) through :

- HACH : HACH LANGE GmbH
- WK : WK Electric Co., Ltd.

Calibrated by : Ailsara Ma

Approved by :

Imron Rattanaylum / Technical Manager

The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



SCaL

Details of Calibration

1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Standard Solution	4.005	C02994	1777	5-Sep-2024
Standard Solution	7.000	C03007	1787	17-Oct-2024
Standard Solution	10.012	C02953	1735	29-Apr-2024
Temperature/Electrical Calibrator	MC2-TE	14987	WK2106-299-223	5-Jun-2024
Digital Thermometer With Sensor	DP-77	1.360896	22SDTH005	8-Aug-2023

2. The results reported in this certificate refer to the condition of the instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.

3. This certificate is not certified any commercial transaction

4. Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration :

1. Electrical Measurement

Applied Voltage (mV)	pH meter Reading		Correction (mV)	Uncertainty (± mV)
	(mV)	(pH)		
177.48	177.5	3.70	-0.02	0.17
0.00	0.0	6.70	0.00	0.13
-177.48	-177.4	9.80	-0.08	0.17

2. Before Sample Test Measurement

Standard Buffer Solutions (pH)	pH meter Reading		Correction (pH)	Uncertainty (± pH)
	(pH)	(mV)		
4.005	3.97	159.2	0.035	0.0090
6.999	6.98	-15.4	0.019	0.013
10.012	9.95	-188.0	0.062	0.036

3. After Sample Test Measurement

Standard Buffer Solutions (pH)	pH meter Reading		Correction (pH)	Uncertainty (± pH)
	(pH)	(mV)		
4.005	3.97	158.9	0.035	0.0090
6.999	7.01	-17.5	-0.011	0.013
10.012	9.98	-187.6	0.032	0.036

4. Temperature Measurement

Cal Point (°C)	Standard Temperature (°C)	UUC Reading (°C)	Correction (°C)	Uncertainty (± °C)
25	25.032	25.1	-0.068	0.11

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k = 2, providing a level of confidence of approximately 95%

... End ...

Certificate No. : 23CH0203

CSR No. : A088/04367

Page : 2 of 2



CARIBRATION CERTIFICATE

Issued Date : 22-Oct-2024

Certificate No. : 24CH0526

CSR No. : A163/08133

Page. : 1 of 2

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Kamchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Chemical Laboratory
Instrument Name : pH meter
Manufacturer : Mettler Toledo
Model : S220
Serial No. : B625631849
ID No. : SGK_CL0030
Electrode No. : 2281592
Received Date : 19-Oct-2024
Calibrated Date : 19-Oct-2024
Ambient Temperature : (25 ± 3) °C
Relative Humidity : (55 ± 15) %

REVIEW BY Ananta B.
APPROVED BY Kanitta H.
NEXT CAL. DATE 19/14/26

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WI.008 based on direct measurement by using certified reference Material (CRM)

The Southern Calibration Service Co.,Ltd.calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and /or national standards which realize the units of measurement according to the International System of Unit (SI) through :

- CPAchem : CPAchem Ltd
- WK : WK Electric Co., Ltd.
- SCAL : Sounthern Calibration Service Co., Ltd.,

Calibrated by : Alisara Ma

Approved by :

Imron Rattanaylum / Technical Manager



The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 24CH0526

CSR No. : A163/08133

Page. : 2 of 2

Details of Calibration

1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Standard Solution	4.000	61310674	1042701	26-Oct-2025
Standard Solution	7.000	61314184	1042700	26-Oct-2025
Standard Solution	10.01	61313804	1042702	26-Oct-2025
Temperature/Electrical Calibrator	MC2-TE	14987	WK2106-299-223	31-May-2025
Digital Thermometer With Sensor	DP-77	1.360896	24SDTH005	7-Aug-2025

2. The results reported in this certificate refer to the condition of the Instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.

3. This certificate is not certified any commercial transaction

4. Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration :

1. Electrical Measurement

Applied Voltage (mV)	pH meter Reading		Correction (mV)	Uncertainty (± mV)
	(mV)	(pH)		
177.48	177.5	4.00	-0.02	0.17
0.00	0.0	7.00	0.00	0.13
-177.48	-177.5	10.00	0.02	0.17

2. Before Sample Test Measurement

Standard Buffer Solutions (pH)	pH meter Reading		Correction (pH)	Uncertainty (± pH)
	(pH)	(mV)		
4.007	3.99	178.1	0.017	0.0092
6.976	7.02	1.7	-0.044	0.019
10.009	9.96	-169.6	0.049	0.038

3. After Sample Test Measurement

Standard Buffer Solutions (pH)	pH meter Reading		Correction (pH)	Uncertainty (± pH)
	(pH)	(mV)		
4.007	3.99	177.7	0.017	0.0092
6.976	7.01	3.7	-0.034	0.019
10.009	10.00	-169.0	0.009	0.038

4. Temperature Measurement

Cal Point (°C)	Standard Temperature (°C)	UUC Reading (°C)	Correction (°C)	Uncertainty (± °C)
25	25.021	25.0	0.02	0.060

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k = 2, providing a level of confidence of approximately 95%

... End ...

Southern Calibration Service Co., Ltd.

669/35 Karnjanavanit Rd., Banpru, Hatyai, Songkla 90250 Thailand
Tel : 08 1599 0417 Fax : 0 7480 5133 Email : s.calibration@gmail.com www.scal-lab.com



CALIBRATION CERTIFICATE

Issued Date : 18-Jan-2024

Certificate No. : 24MA0199

CSR No. : A123/06123

Page. : 1 of 3

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Kamchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Chemical Laboratory
Instrument Name : Electronic Balance
Manufacturer : Sartorius
Model : MSE224S-100-DU
Serial No. : 34705158
ID No. : SGK_CL0045
Resolution : 0.0001 g
Received Date : 15-Jan-2024
Calibrated Date : 15-Jan-2024
Ambient Temperature : (30 ± 10) °C
Relative Humidity : (50 ± 20) %

REVIEW BY Ananta B
APPROVED BY Kamilla H.
NEXT CAL. DATE 15/11/25

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WI.001 based on UKAS LAB 14 : 2015

The Southern Calibration Service Co., Ltd. calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and/or national standards which realize the units of measurement according to the International System of Unit (SI) through :

- SCaL : Southern Calibration Service Co., Ltd.,

Calibrated by : Hadbordee Dettawee

Approved by :

Imron Rattanaylum / Technical Manager



The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 24MA0199

CSR No. : A123/06123

Page. : 2 of 3

Details of Calibration

1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Standard Weight Set	2 mg - 1 kg	11119514/01	23SWS001	4-Jul-2024

2. The results reported in this certificate refer to the condition of the instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.

3. This certificate is not certified any commercial transaction

4. Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration : (✓) Without Adjustment () After Adjustment

1. Repeatability

Nominal Value (g)	Standard Deviation (g)
20	0.00000
200	0.00000

2. Effect of tare

Nominal Value (g)	Standard Value (g)	Balance Reading (g)	Correction (g)
20	20.0000	20.0000	0.0000
40	40.0001	40.0000	0.0001
60	60.0000	60.0001	-0.0001
80	80.0001	80.0001	-0.0001
100	100.0000	100.0000	0.0000



Certificate No. : 24MA0199

CSR No. : A123/06123

Page. : 3 of 3


Result of Calibration :

3. Off-centre loading

A mass approximately 100g was placed on a pan and moved to various position .

The balance reading obtained are given in the table.

Position					Maximum Difference (g)
1	2	3	4	5	
100.0000	100.0001	100.0000	100.0000	100.0000	0.0001

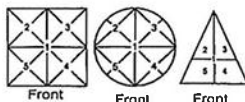

4. Departure from nominal value

Nominal Value (g)	Standard Value (g)	UUC Reading (g)	Correction (g)	Uncertainty ($\pm$ g)	Coverage Factor (k)
0	0.0000	0.0000	0.0000	0.00008	2.0
0.01	0.0100	0.0100	0.0000	0.00008	2.0
0.1	0.1000	0.1000	0.0000	0.00008	2.0
0.5	0.5000	0.5000	0.0000	0.00008	2.0
1	1.0000	1.0000	0.0000	0.00008	2.0
2	2.0000	2.0000	0.0000	0.00008	2.0
5	5.0000	5.0000	0.0000	0.00009	2.0
10	10.0000	10.0000	0.0000	0.00009	2.0
20	20.0000	20.0000	0.0000	0.00009	2.0
50	50.0000	50.0000	0.0000	0.00011	2.0
100	100.0000	100.0000	0.0000	0.00016	2.0
120	120.0000	120.0000	0.0000	0.00024	2.0
140	140.0001	140.0000	0.0001	0.00024	2.0
160	160.0000	160.0000	0.0000	0.00026	2.0
180	180.0000	180.0000	0.0000	0.00029	2.0
200	200.0000	200.0000	0.0000	0.00030	2.0

- UUC = Unit Under Calibration

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95%

...End...



Southern Calibration Service Co., Ltd.

669/35 Karnjanavanit Rd., Banpru, Hatyai, Songkla 90250 Thailand
Tel : 08 1599 0417 Fax : 0 7480 5133 Email : s.calibration@gmail.com www.scal-lab.com



CALIBRATION CERTIFICATE

Issued Date : 1-May-2023

Certificate No. : 23TH1728

CSR No. : A088/04367

Page. : 1 of 3

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Karnchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Chemical Laboratory
Instrument Name : Hot Air Oven
Manufacturer : Memmert
Model : UF110
Serial No. : B416.3387
ID No. : SGK_CL0024
Resolution : 0.1 °C
Received Date : 28-Apr-2023
Calibrated Date : 28-Apr-2023
Ambient Temperature : (30 $\pm$ 10) °C
Relative Humidity : (50 $\pm$ 30) %

REVIEW BY Ananta B.
APPROVED BY Kamrta H.
NEXT CAL. DATE 28/11/2024

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WI.012 based on GLA - 20

The Southern Calibration Service Co.,Ltd.calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and/or national standards which realize the units of measurement

according to the International System of Unit (SI) through :

- SCaL : Southern Calibration Service Co., Ltd.,

Calibrated by : Ibrahim Saleemin

Approved by :

Imron Rattanaylum / Technical Manager

The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 23TH1728
CSR No. : A088/04367
Page. : 2 of 3

Details of Calibration


1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Data Acquisition/Switch Unit	34970A	MY58009813	22SDAT004	24-May-2023

2. The results reported in this certificate refer to the condition of the instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.
3. This certificate is not certified any commercial transaction
4. Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration : (✓) Without Adjustment () After Adjustment

1. Sensor Installation Diagram

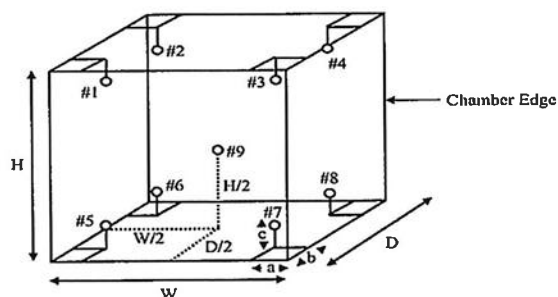

Sensor Installation Details

a = 5.0 cm
b = 5.0 cm
c = 5.0 cm

Dimension of the chamber

W = 40.0 cm
H = 40.0 cm
D = 33.0 cm



Certificate No. : 23TH1728
CSR No. : A088/04367
Page. : 3 of 3

Result of Calibration :

2. Temperature Measurement Accuracy Test

The measurement results of the Hot Air Oven and associates are reported in the manner as shown below

Cal point (°C)	Measured Standard Temperature At Spread Locations (°C)									Uncertainty (± °C)
	#1	#2	#3	#4	#5	#6	#7	#8	Ref. 9	
40	40.48	40.28	40.28	39.91	40.17	40.09	39.93	40.27	39.89	0.36
70	70.36	70.23	70.58	69.74	69.99	69.92	69.86	70.13	70.04	0.36
103	103.19	103.12	103.46	103.37	103.10	103.54	103.43	103.06	103.40	0.36
104	104.31	104.23	104.62	103.77	104.12	104.06	103.90	104.20	104.56	0.36
105	105.07	105.03	105.48	105.27	105.12	105.01	105.01	105.00	104.96	0.36
180	180.31	180.00	180.00	180.07	180.18	180.05	180.01	180.10	180.24	0.41

3. Performance Result

The performance of the Hot Air Oven are reported as shown below

Cal point (°C)	UUC Setting (°C)	UUC Reading (°C)	Temperature Stability (± °C)	Temperature Uniformity (°C)	Overall Variation (°C)
40	40.0	40.0	0.20	0.70	0.72
70	70.0	70.0	0.20	0.60	0.94
103	103.0	103.0	0.20	0.43	0.54
104	104.0	104.0	0.10	0.79	0.88
105	105.0	105.0	0.10	0.59	0.69
180	180.0	180.0	0.10	0.38	0.38

- UUC = Unit Under Calibration

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%

... End ...



CARIBRATION CERTIFICATE

Issued Date : 22-Oct-2024

Certificate No. : 24TH4295

CSR No. : A163/08133

Page. : 1 of 3

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Kamchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Chemical Laboratory

Instrument Name : Hot Air Oven

Manufacturer : Memmert

Model : UF110

Serial No. : B416.3387

ID No. : SGK_CL0024

Resolution : 0.1 °C

Received Date : 19-Oct-2024

Calibrated Date : 19-Oct-2024

Ambient Temperature : (30 ± 10) °C

Relative Humidity : (50 ± 30) %

REVIEW BY Ananta B.
APPROVED BY Kanitta H.
NEXT CAL. DATE 19/04/2026

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WI.012 based on GLA - 20

The Southern Calibration Service Co.,Ltd.calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and /or national standards which realize the units of measurement
according to the International System of Unit (SI) through :

- TISTR : Thailand Institute of Scientific and Technological Research

Calibrated by : Ibrorhim Saleemin

Approved by :

Imron Rattianaylum / Technical Manager



The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 24TH4295

CSR No. : A163/08133

Page. : 2 of 3

Details of Calibration


1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Data Acquisition/Switch Unit	34970A	MY58009813	PSL-T0707-1/67	22-May-2025

- The results reported in this certificate refer to the condition of the instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.
- This certificate is not certified any commercial transaction
- Condition of item : normal condition , no indication for any damage or malfunction

Result of Calibration : (✓) Without Adjustment () After Adjustment

1. Sensor Installation Diagram

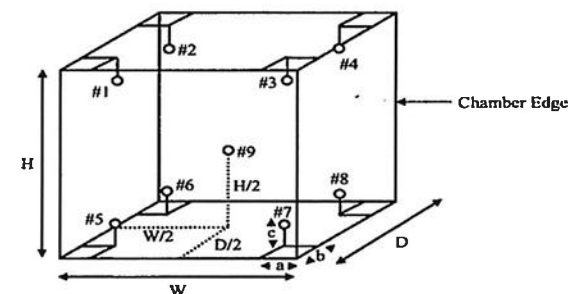

Sensor Installation Details

a = 5.0 cm
b = 5.0 cm
c = 5.0 cm

Dimension of the chamber

W = 55.0 cm
H = 48.0 cm
D = 40.0 cm



Certificate No. : 24TH4295

CSR No. : A163/08133

Page. : 3 of 3

Result of Calibration :**2. Temperature Measurement Accuracy Test**

The measurement results of the Hot Air Oven and associates are reported in the manner as shown below

Cal point (°C)	Measured Standard Temperature At Spread Locations (°C)									Uncertainty (±°C)
	#1	#2	#3	#4	#5	#6	#7	#8	Ref. 9	
40	40.36	40.40	40.51	40.43	40.05	40.24	40.09	40.14	39.75	0.38
70	70.27	70.30	70.45	70.24	70.24	70.43	70.29	70.30	69.95	0.36
103	102.94	102.90	103.55	102.96	103.22	103.14	103.10	103.01	102.88	0.36
104	104.15	103.99	104.27	104.06	104.09	104.23	104.26	104.15	103.90	0.36
105	105.04	104.90	105.05	104.87	104.91	104.80	104.82	104.98	104.70	0.36
180	179.19	178.93	179.82	179.10	179.27	179.68	179.12	179.73	179.12	0.41

3. Performance Result

The performance of the Hot Air Oven are reported as shown below

Cal point (°C)	UUC Setting (°C)	UUC Reading (°C)	Temperature Stability (±°C)	Temperature Uniformity (°C)	Overall Variation (°C)
40	40.0	40.0	0.20	0.84	0.84
70	70.0	70.0	0.10	0.59	0.59
103	103.0	103.0	0.20	0.73	0.74
104	104.0	104.0	0.20	0.47	0.56
105	105.0	105.0	0.20	0.44	0.46
180	180.0	180.0	0.50	0.86	1.11

- UUC = Unit Under Calibration

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$,
providing a level of confidence of approximately 95%

... End ...

**Internal Calibration**

Equipment : Turbidity Meter Manufacture : HACH
ID No. : SGK FS0045 Model : 2100Q/QIS
Calibrate Date : May 22, 2024 Serial No. : 19010C073443

Calibration Point	1st (NTU)	2nd (NTU)	3rd (NTU)	AVG (NTU)	Specifications	Evaluate
Standard 20 NTU	19.8	20.1	20.0	20.0	19 to 21 NTU	Pass
Standard 100 NTU	99.5	99.3	100	99.6	95 to 105 NTU	Pass
Standard 800 NTU	798	799	798	798.3	760 to 840 NTU	Pass

REVIEW BY	Somsak J.
APPROVED BY	Kanida H.
NEXT CAL. DATE	22 May 25

Calibrated by Somsak J.
(Scientist 2)

Approved by : Kanida H.
(Section Head)

HACH COMPANY



An ISO 9001 Certified Company

P.O.Box 389
Loveland, CO 80539
(970) 669-3050

Certificate of Analysis

This is a Component of 2659405-TH lot A4025

Page 1

COMMODITY:STABLCAL|sup|TS|sup0 FORMAZIN STANDARD 20 NTU
COMMODITY NUMBER: 2684801 MANUFACTURE DATE: DATE OF ANALYSIS:
LOT NUMBER: A4025 2/6/2024 2/12/2024

TEST	SPECIFICATIONS	RESULTS
Turbidity	19 to 21 NTU	20.5 NTU

The expiration date is May 2025

Formazin and StablCal® solutions provided by Hach are not NIST traceable because the NIS does not carry turbidity standards. However, the use of Formazin and StablCal® as use in Hach method 8195 are accepted by the EPA as a primary standard to be used in the calibration of turbidity instruments.

Certified by

Scott Als
Analytical Services Chemist

HACH COMPANY



An ISO 9001 Certified Company

P.O.Box 389
Loveland, CO 80539
(970) 669-3050

Certificate of Analysis

This is a Component of 2659405-TH lot A4025

Page 1

COMMODITY:STABLCAL|sup|TS|sup0 FORMAZIN STANDARD 100 NTU
COMMODITY NUMBER: 2684901 MANUFACTURE DATE: DATE OF ANALYSIS:
LOT NUMBER: A4018 1/29/2024 2/2/2024

TEST	SPECIFICATIONS	RESULTS
Turbidity	95 to 105 NTU	102.0 NTU

The expiration date is Apr 2025

Formazin and StablCal® solutions provided by Hach are not NIST traceable because the NIS does not carry turbidity standards. However, the use of Formazin and StablCal® as use in Hach method 8195 are accepted by the EPA as a primary standard to be used in the calibration of turbidity instruments.

Certified by

Scott Als
Analytical Services Chemist

HACH COMPANY



An ISO 9001 Certified Company

P.O.Box 389
Loveland, CO 80539
(970) 669-3050

Certificate of Analysis

This is a Component of 2659405-TH lot A4025

Page 1
COMMODITY: STABLCAL|sup|TS|sup0 FORMAZIN STANDARD 800 NTU
COMMODITY NUMBER: 2660501 MANUFACTURE DATE: DATE OF ANALYSIS:
LOT NUMBER: A4017 2/1/2024 2/9/2024

TEST	SPECIFICATIONS	RESULTS
Turbidity	760 to 840 NTU	818.0 NTU

The expiration date is May 2025

Formazin and StablCal® solutions provided by Hach are not NIST traceable because the NIS does not carry turbidity standards. However, the use of Formazin and StablCal® as used in Hach method 8195 are accepted by the EPA as a primary standard to be used in the calibration of turbidity instruments.

Certified by

Scott Als
Analytical Services Chemist



Southern Calibration Service Co., Ltd.

669/35 Karnjanavanit Rd., Banpru, Hatyai, Songkla 90250 Thailand
Tel : 08 1599 0417 Fax : 0 7480 5133 Email : s.calibration@gmail.com www.seul-lab.com



CALIBRATION CERTIFICATE

Issued Date : 16-Jul-2023

Certificate No. : 23TH3096

CSR No. : A095/04743

Page. : 1 of 3

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Karnchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Chemical Laboratory
Instrument Name : Incubator
Manufacturer : Memmert
Model : ICP750
Serial No. : F816.0063
ID No. : SGK_CL0028
Resolution : 0.1 °C
Received Date : 13-Jul-2023
Calibrated Date : 13-Jul-2023
Ambient Temperature : (30 ± 10) °C
Relative Humidity : (50 ± 30) %

REVIEW BY Ananta B.
APPROVED BY Kanytha H.
NEXT CAL. DATE 13/01/25

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WI.012 based on GLA - 20

The Southern Calibration Service Co., Ltd. calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and /or national standards which realize the units of measurement according to the International System of Unit (SI) through :

- SCAL : Southern Calibration Service Co., Ltd.,

Calibrated by : Ibrahim Saleemin

Approved by :

Imron Rattanaylum / Technical Manager

The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 23TH3096
CSR No. : A095/04743
Page. : 2 of 3

Details of Calibration


1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Data Acquisition/Switch Unit	34970A	MY58009813	23SDAT004	23-May-2024

2. The results reported in this certificate refer to the condition of the instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.
3. This certificate is not certified any commercial transaction
4. Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration : (✓) Without Adjustment () After Adjustment

1. Sensor Installation Diagram

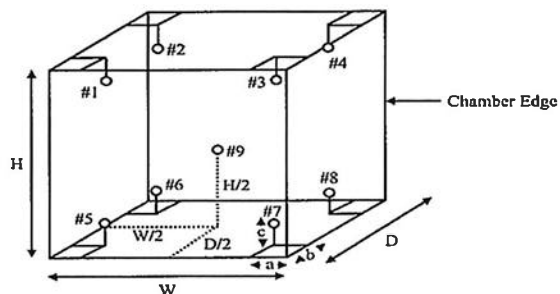

Sensor Installation Details

a = 5.0 cm
b = 5.0 cm
c = 5.0 cm

Dimension of the chamber

W = 40.0 cm
H = 40.0 cm
D = 33.0 cm



Certificate No. : 23TH3096
CSR No. : A095/04743
Page. : 3 of 3

Result of Calibration :

2. Temperature Measurement Accuracy Test

The measurement results of the incubator and associates are reported in the manner as shown below

Cal point (°C)	Measured Standard Temperature At Spread Locations (°C)									Uncertainty (± °C)
	#1	#2	#3	#4	#5	#6	#7	#8	Ref. 9	
20	20.10	20.04	20.03	19.97	20.08	20.23	20.10	19.94	20.07	0.38

3. Performance Result

The performance of the Incubator are reported as shown below

Cal point (°C)	UUC Setting (°C)	UUC Reading (°C)	Temperature Stability (± °C)	Temperature Uniformity (°C)	Overall Variation (°C)
20	20.0	20.0	0.14	0.17	0.32

- UUC = Unit Under Calibration

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%

... End ...



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3 : EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000 FAX. 0-2719-9484

Certificate of Testing

Cert.No.: 24TW96
Page.: 1 of 2

Equipment : DO Meter
Manufacturer : YSI
Model : 5000
Serial No. : 17B101473
ID No. : SGK_CL0073
Received Date : 17 May 2024
Test Date : 21 May 2024
Reference : 2405-0608DSC-1
Submitted by : ALS Laboratory Group (Thailand) Co.,Ltd.
Songkhla Branch.
114/1 Moo 8 Kamchanawanich Rd., T.Ban Phru,
A.Hat Yai, Songkhla 90250 Thailand

Laboratory Condition : Temperature (25 ± 5) °C
Humidity (50 ± 20) %
Test Procedure : In - house method : CP-CH9
by Comparison Technique with Azide Modification Method

Tested by : Walalak Sirithean

Approved by : Saithip
Approved Signatory

() Unnopphol Harachai
() Ponpan Paipim
(✓) Saithip Meangmai

Issue Date : 21 May 2024



Cert.No.: 24TW96
Page.: 2 of 2

Condition of this result of calibration

1. Reference Standard Instruments :

This certification is traceable to the International System of Unit through the reference standards laboratory of Industrial Calibration Center, Technology Promotion Association (Thailand-Japan).

Instruments	Serial No.	ID No.	Certificate No.	Due Date
1. Burette	-	130BU10	23CG1172	22 Mar 2025
2. Balance	14233821	110RC001	23MM405	16 July 2024

2. Standard Material :-

Material	Manufacturer	Lot.No.	Assay
Sodium Thiosulfate pentahydrate	Merck	AM1763316	100.2%

Result : Dissolved Oxygen Meter Adjustment With Air 100 %
Dissolved Oxygen Probe No.: 17B100103

Titration Method (Azide Modification Method) (mg/L)	DO Meter Reading (mg/L)	Standard Deviation (mg/L)
8.18	8.18	0.0071

This report was certified only for the instrument we tested.It is allowable to use for study
Intend to use for advertising and referral purpose is prohibited.This report may not be reproduced
other in full,without written approval of the laboratory

-o0o-




TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert. No.: 24LM77
Page.: 1 of 2

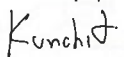Equipment : DO Meter with Sensor
Manufacturer : YSI
Model : 5000-115
Serial No. : 17B101473
ID No. : SGK_CL0073
Submitted by : ALS Laboratory Group (Thailand) Co.,Ltd.
Songkhla Branch.
114/1 Moo 8 Kamchanawanich Rd.,
T.Ban Phru, A.Hat Yai,
Songkhla 90250 Thailand
Location : TPA On Site Calibration Laboratory
Received Order : 17 May 2024
Calibrated Date : 27 May 2024
Ambient Temperature : (26 ± 10) °C
Relative Humidity : (50 ± 30) %
AC Line Voltage : (220 ± 22) V
Calibrated by : Khit Ruttanaprapachai
Approved by : 
Approved Signatory
() Ponpan Paipim
() Suwit Imjai
(✓) Kunchit Promprat
Issue Date : 28 May 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Equipment : DO Meter with Sensor
Condition As-Received : Used Item
Reference : 2405-0608DSC-2
Procedure Used :-

Cert. No.: 24LM77
Page.: 2 of 2

Calibration were conducted using in-house calibration procedure CP-OT01 according to comparison with Industrial Platinum Resistance Thermometer (IPRT) into Temperature Bath.
The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Digital Thermometer	2188080	2311216	TPA	11 Oct 2024

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function : Temperature measurement.

This instrument was connected with temperature sensor, S/N.: 17B100103

Calibration Point (°C)	Immersion Depth (mm)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty (± °C)	Coverage Factor k
20.00	60	20.005	19.79	-0.215	0.15	2.00

UUC* : Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-



Southern Calibration Service Co., Ltd.

669/35 Karnjanavanit Rd., Banpru, Hatyai, Songkla 90250 Thailand
Tel : 08 1599 0417 Fax : 0 7480 5133 Email : s.calibration@gmail.com www.scal-lab.com



CALIBRATION CERTIFICATE

Issued Date : 16-Jul-2023

Certificate No. : 23TH3097

CSR No. : A095/04743

Page : 1 of 3

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Karnchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Chemical Laboratory
Instrument Name : Water Bath
Manufacturer : Memmert
Model : WNE29
Serial No. : L616.0538
ID No. : SGK_CL0035
Resolution : 0.1 °C
Received Date : 13-Jul-2023
Calibrated Date : 13-Jul-2023
Ambient Temperature : (30 ± 10) °C
Relative Humidity : (50 ± 30) %

REVIEW BY Ananta B.
APPROVED BY Kanitha H.
NEXT CAL. DATE 13/01/25

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WI.014 based on ASTM E 715 : 1980
(reapproved 2001)

The Southern Calibration Service Co., Ltd. calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and/or national standards which realize the units of measurement
according to the International System of Unit (SI) through :

- SCal : Southern Calibration Service Co., Ltd.,

Calibrated by : Ibrahim Saleemin

Approved by :

Imron Rattanaylum / Technical Manager

Certificate No. : 23TH3097

CSR No. : A095/04743

Page : 2 of 3



Details of Calibration

1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Data Acquisition/Switch Unit	34970A	MY58009813	23SDAT004	23-May-2024


2. The results reported in this certificate refer to the condition of the instrument on the date of calibration
and carry no implication regarding the longterm stability of instrument.

3. This certificate is not certified any commercial transaction

4. Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration : (✓) Without Adjustment () After Adjustment

1. Sensor Installation Diagram

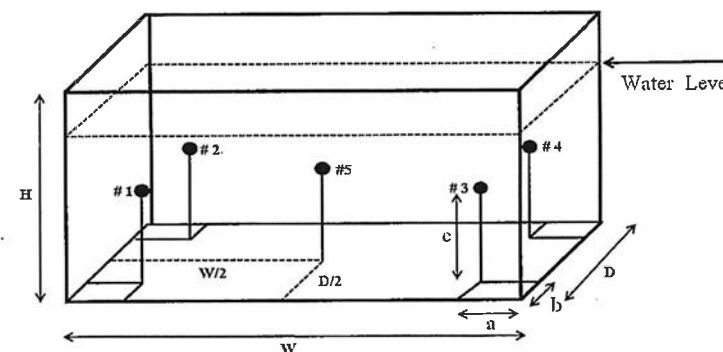

Sensor Installation Details

a = 5 cm
b = 5 cm
c = 5 cm

Dimension of the chamber

W = 45 cm
H = 30 cm
D = 35 cm

The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 23TH3097
CSR No. : A095/04743
Page. : 3 of 3

Result of Calibration :

2. Temperature Measurement Accuracy Test

The measurement results of the Water Bath and associates are reported in the manner as shown below

Cal point (°C)	Measured Standard Temperature At Spread Locations (°C)					Uncertainty (± °C)
	#1	#2	#3	#4	Ref.5	
80	79.17	79.47	79.43	79.25	79.38	0.14

3. Performance Result


The performance of the Water Bath are reported as shown below

Cal point (°C)	UUC Setting (°C)	UUC Reading (°C)	Temperature Stability (± °C)	Temperature Uniformity (°C)	Overall Variation (°C)
80	80.0	80.0	0.24	0.38	0.38

- UUC = Unit Under Calibration

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%

... End ...

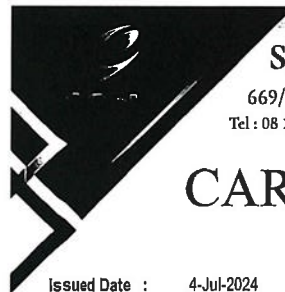

Southern Calibration Service Co., Ltd.
669/35 Karnjanavanit Rd., Banpru, Hatyai, Songkla 90250 Thailand
Tel : 08 1599 0417 Fax : 0 7480 5133 Email : s.calibration@gmail.com www.scal-lab.com

CARIBRATION CERTIFICATE



NSC-TISI-TIS17025
CALIBRATION 0254

Issued Date : 4-Jul-2024

Certificate No. : 24TH2760
CSR No. : A150/07473
Page. : 1 of 3

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Kamchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Microbiological Laboratory
Instrument Name : Autoclave
Manufacturer : TOMY
Model : SX-700
Serial No. : 52134079
ID No. : SGK_ML0001
Resolution : 1 °C
Received Date : 1-Jul-2024
Calibrated Date : 1-Jul-2024
Ambient Temperature : (30 ± 10) °C
Relative Humidity : (50 ± 30) %

REVIEW BY Wichapol P.
APPROVED BY Kanida H.
NEXT CAL. DATE 1/01/26

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WI.013 based on BS 2646 : 1993 (part 5)
The Southern Calibration Service Co.,Ltd.calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and /or national standards which realize the units of measurement according to the International System of Unit (SI) through :

- SCaL : Sounthern Calibration Service Co., Ltd.,

Calibrated by : Ibrorhim Saleemin

Approved by :

Imron Rattanaylum / Technical Manager



The uncertainties are for a confidence probability of approximately 95%
This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 24TH2760

CSR No. : A150/07473

Page. : 2 of 3

Details of Calibration


1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Data logger With Sensor	GL240	C90432223	24SDAT005	2-May-2025

- The results reported in this certificate refer to the condition of the instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.
- This certificate is not certified any commercial transaction
- Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration.: (✓) Without Adjustment () After Adjustment

1. Sensor Installation Diagram

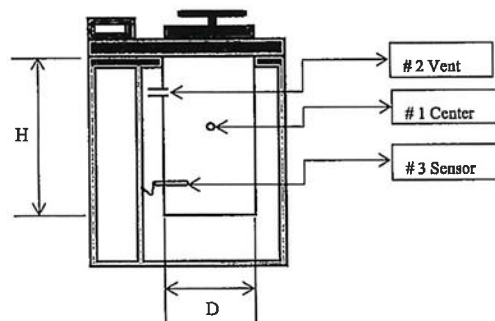

Chamber Diameter (D) : 30 cm

Chamber Height (H) : 70 cm



Certificate No. : 24TH2760

CSR No. : A150/07473

Page. : 3 of 3

Result of Calibration.:

2. Temperature Measurement Accuracy Test

The measurement results of the Autoclave and associates are reported in the manner as shown below

Cal point (°C)	Measured Standard Temperature At Spread Locations (°C)			Pressure Reading	Uncertainty (±°C)
	Center #1	Vent #2	Sensor #3		
115	116.3	116.4	116.3	0.07 MPa	0.76
118	119.2	119.2	119.2	0.09 MPa	0.76
121	121.8	121.5	121.6	0.11 MPa	0.76

3. Performance Result

The performance of the Autoclave are reported as shown below

Cal point (°C)	UUC Setting (°C)	UUC Reading (°C)	Temperature Stability (±°C)	Temperature Uniformity (°C)	Overall Variation (°C)
115	115	115	0.10	0.17	0.17
118	118	118	0.10	0.10	0.10
121	121	121	0.70	0.50	0.84

- UUC = Unit Under Calibration

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%

... End ...



CARIBRATION CERTIFICATE

Issued Date : 18-Jan-2024

Certificate No. : 24TH0186

CSR No. : A123/06123

Page. : 1 of 3

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Karnchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Microbiological Laboratory
Instrument Name : Incubator
Manufacturer : Memmert
Model : ICP750
Serial No. : F816.0061
ID No. : SGK_ML0013
Resolution : 0.1 °C
Received Date : 15-Jan-2024
Calibrated Date : 15-Jan-2024
Ambient Temperature : (30 ± 10) °C
Relative Humidity : (50 ± 30) %

REVIEW BY Wichapol R.
APPROVED BY Kamilla H.
NEXT CAL. DATE 15/07/25

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WI.012 based on GLA - 20

The Southern Calibration Service Co., Ltd. calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and /or national standards which realize the units of measurement
according to the International System of Unit (SI) through :

- SCS : Southern Calibration Service Co., Ltd.,

Calibrated by : Ibrorhim Saleemin

Approved by :

Imron Rattanaylum / Technical Manager



The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 24TH0186

CSR No. : A123/06123

Page. : 2 of 3

Details of Calibration

1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Data Acquisition/Switch Unit	34970A	MY58009813	23SDAT004	23-May-2024


2. The results reported in this certificate refer to the condition of the instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.
3. This certificate is not certified any commercial transaction
4. Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration :

(✓) Without Adjustment

() After Adjustment

1. Sensor Installation Diagram

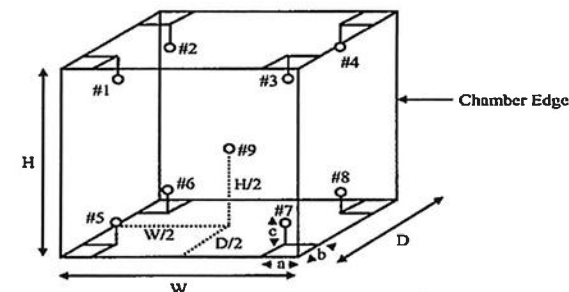

Sensor Installation Details

a = 5.0 cm
b = 5.0 cm
c = 5.0 cm

Dimension of the chamber

W = 104.0 cm
H = 120.0 cm
D = 60.0 cm



Certificate No. : 24TH0186
CSR No. : A123/06123
Page. : 3 of 3

Result of Calibration :

2. Temperature Measurement Accuracy Test

The measurement results of the Incubator and associates are reported in the manner as shown below

Cal point (°C)	Measured Standard Temperature At Spread Locations (°C)									Uncertainty (±°C)
	#1	#2	#3	#4	#5	#6	#7	#8	Ref. 9	
35	35.04	34.97	34.90	34.89	35.05	34.91	34.99	34.92	34.94	0.38

3. Performance Result

The performance of the Incubator are reported as shown below

Cal point (°C)	UUC Setting (°C)	UUC Reading (°C)	Temperature Stability (±°C)	Temperature Uniformity (°C)	Overall Variation (°C)
35	35.0	35.0	0.10	0.23	0.25

- UUC = Unit Under Calibration

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%

... End ...

Southern Calibration Service Co., Ltd.

669/35 Karnjanavanit Rd., Banpru, Hatyai, Songkla 90250 Thailand
Tel : 08 1599 0417 Fax : 0 7480 5133 Email : s.calibration@gmail.com www.scal-lab.com



CALIBRATION CERTIFICATE

Issued Date : 5-Jul-2024

Certificate No. : 24CH0342

CSR No. : A150/07473

Page. : 1 of 2

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Karnchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Microbiological Laboratory
Instrument Name : pH meter
Manufacturer : Sartorius
Model : PB-10
Serial No. : C07160695
ID No. : SGK_ML0016
Electrode No. : P20067003
Received Date : 2-Jul-2024
Calibrated Date : 2-Jul-2024
Ambient Temperature : (25 ± 3) °C
Relative Humidity : (55 ± 15) %

REVIEW BY	Wichapol D.
APPROVED BY	Kanika H.
NEXT CAL. DATE	2/01/26

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WI.008 based on direct measurement by using certified reference Material (CRM)

The Southern Calibration Service Co., Ltd. calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and/or national standards which realize the units of measurement according to the International System of Unit (SI) through :

- CPAchem : CPAchem Ltd
- WK : WK Electric Co., Ltd.

Calibrated by : Alisara Ma

Approved by :

Imron Rattanaylum / Technical Manager



The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 24CH0342
CSR No. : A150/07473
Page : 2 of 2

Details of Calibration

1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Standard Solution	4.000	61278486	986239	10-May-2025
Standard Solution	7.000	61281486	986241	10-May-2025
Standard Solution	10.01	61287532	986240	10-May-2025
Temperature/Electrical Calibrator	MC2-TE	14987	WK2106-299-223	31-May-2025

2. The results reported in this certificate refer to the condition of the instrument on the date of calibration

and carry no implication regarding the longterm stability of instrument.

3. This certificate is not certified any commercial transaction

4. Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration :

1. Electrical Measurement

Applied Voltage (mV)	pH meter Reading		Correction (mV)	Uncertainty (± mV)
	(mV)	(pH)		
177.48	177.5	3.80	-0.02	0.17
0.00	0.0	6.90	0.00	0.13
-177.48	-177.5	10.00	0.02	0.17

2. Before Sample Test Measurement

Standard Buffer Solutions (pH)	pH meter Reading		Correction (pH)	Uncertainty (± pH)
	(pH)	(mV)		
4.008	4.07	164.0	-0.062	0.011
6.985	7.02	-7.1	-0.035	0.020
10.010	10.03	-182.1	-0.020	0.037

3. After Sample Test Measurement

Standard Buffer Solutions (pH)	pH meter Reading		Correction (pH)	Uncertainty (± pH)
	(pH)	(mV)		
4.008	4.02	164.2	-0.012	0.011
6.985	7.00	-0.9	-0.015	0.020
10.010	10.01	-170.7	0.000	0.037

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%

... End ...



Southern Calibration Service Co., Ltd.

669/35 Karnjanavanit Rd., Banpru, Hatyai, Songkla 90250 Thailand
Tel : 08 1599 0417 Fax : 0 7480 5133 Email : s.calibration@gmail.com www.scal-lab.com



CALIBRATION CERTIFICATE

Issued Date : 4-Jul-2024

Certificate No. : 24TH2782

CSR No. : A150/07474

Page : 1 of 3

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Kamchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Microbiological Laboratory
Instrument Name : Liquid Bath
Manufacturer : Memmert
Model : WPE45
Serial No. : L716.0558
ID No. : SGK_ML0021
Resolution : 0.1 °C
Received Date : 1-Jul-2024
Calibrated Date : 1-Jul-2024
Ambient Temperature : (30 ± 10) °C
Relative Humidity : (50 ± 30) %

REVIEW BY Wichapol R.
APPROVED BY Kanitta H.
NEXT CAL. DATE 1/01/26

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WI.014 based on ASTM E 715 : 1980 (reapproved 2001)

The Southern Calibration Service Co., Ltd. calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and /or national standards which realize the units of measurement according to the International System of Unit (SI) through :

- TISTR : Thailand Institute of Scientific and Technological Research

Calibrated by : Ibrorhim Saleemin

Approved by :

Imron Rattanaylum / Technical Manager



The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 24TH2782
CSR No. : A150/07474
Page. : 2 of 3

Details of Calibration


1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Data Acquisition/Switch Unit	34970A	MY58009813	PSL-T0707-1/67	22-May-2025

2. The results reported in this certificate refer to the condition of the Instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.
3. This certificate is not certified any commercial transaction
4. Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration : (✓) Without Adjustment () After Adjustment

1. Sensor Installation Diagram

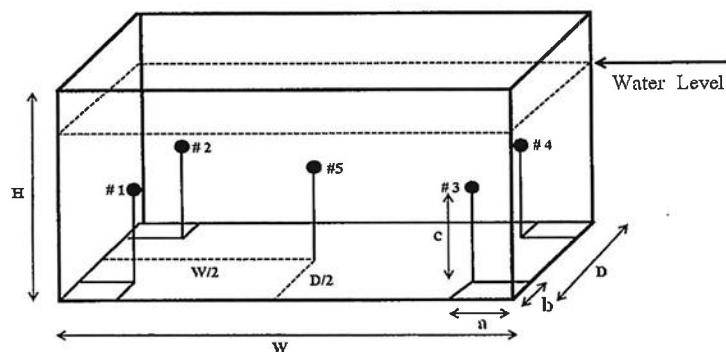

Sensor Installation Details

a = 5 cm
b = 5 cm
c = 5 cm

Dimension of the chamber

W = 45 cm
H = 30 cm
D = 35 cm



Certificate No. : 24TH2782
CSR No. : A150/07474
Page. : 3 of 3

Result of Calibration :

2. Temperature Measurement Accuracy Test

The measurement results of the Liquid Bath and associates are reported in the manner as shown below

Cal point (°C)	Measured Standard Temperature At Spread Locations (°C)					Uncertainty (± °C)
	#1	#2	#3	#4	Ref.5	
44.5	44.48	44.46	44.51	44.49	44.47	0.14

3. Performance Result

The performance of the Liquid Bath are reported as shown below

Cal point (°C)	UUC Setting (°C)	UUC Reading (°C)	Temperature Stability (± °C)	Temperature Uniformity (°C)	Overall Variation (°C)
44.5	44.5	44.5	0.10	0.14	0.25

- UUC = Unit Under Calibration

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%

... End ...

Southern Calibration Service Co., Ltd.

669/35 Karnjanavanit Rd., Banpru, Hatyai, Songkla 90250 Thailand
Tel : 08 1599 0417 Fax : 0 7480 5133 Email : s.calibration@gmail.com www.scal-lab.com



CALIBRATION CERTIFICATE

Issued Date : 4-Jul-2024

Certificate No. : 24TH2761

CSR No. : A150/07473

Page. : 1 of 3

Customer : ALS Laboratory Group (Thailand) Co., Ltd.
114/1 Moo 8, Kamchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Microbiological Laboratory
Instrument Name : Autoclave
Manufacturer : TOMY
Model : SX-700
Serial No. : 52134080
ID No. : SGK_ML0002
Resolution : 1 °C
Received Date : 1-Jul-2024
Calibrated Date : 1-Jul-2024
Ambient Temperature : (30 ± 10) °C
Relative Humidity : (50 ± 30) %

REVIEW BY Wichapol P.
APPROVED BY Kanitta H.
NEXT CAL. DATE 1/01/26

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WI.013 based on BS 2646 : 1993 (part 5)

The Southern Calibration Service Co., Ltd. calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and /or national standards which realize the units of measurement according to the International System of Unit (SI) through :

- SCaL : Southern Calibration Service Co., Ltd.,

Calibrated by : Ibromhim Saleemin

Approved by :

Imron Rattanaylum / Technical Manager



The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 24TH2761

CSR No. : A150/07473

Page. : 2 of 3

Details of Calibration

1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Data logger With Sensor	GL240	C90432223	24SDAT005	2-May-2025


2. The results reported in this certificate refer to the condition of the instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.

3. This certificate is not certified any commercial transaction

4. Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration : (✓) Without Adjustment () After Adjustment

1. Sensor Installation Diagram

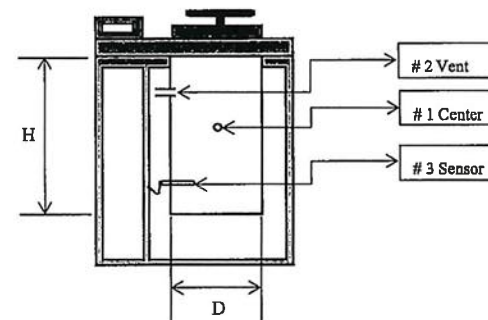

Chamber Diameter (D) : 30 cm

Chamber Height (H) : 70 cm



Certificate No. : 24TH2761

CSR No. : A150/07473

Page. : 3 of 3

Result of Calibration :

2. Temperature Measurement Accuracy Test

The measurement results of the Autoclave and associates are reported in the manner as shown below

Cal point (°C)	Measured Standard Temperature At Spread Locations (°C)			Pressure Reading	Uncertainty (±°C)
	Center #1	Vent #2	Sensor #3		
115	116.3	115.9	116.2	0.07 MPa	0.76
118	119.6	119.4	119.2	0.09 MPa	0.76
121	122.0	122.0	122.1	0.11 MPa	0.76

3. Performance Result

The performance of the Autoclave are reported as shown below

Cal point (°C)	UUC Setting (°C)	UUC Reading (°C)	Temperature Stability (±°C)	Temperature Uniformity (°C)	Overall Variation (°C)
115	115	115	0.10	0.33	0.43
118	118	118	0.30	0.50	0.50
121	121	121	0.20	0.14	0.34

- Operating Time = 900.00 sec

- UUC = Unit Under Calibration

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%

... End ...



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert.No.: 24CG952

Page.: 1 of 2

Equipment : Burette
Capacity : 50 mL
Serial No. : -
ID. No. : BKK_EN0171
Manufacturer : Witeg
Made in : Germany
Submitted by :

ALS Laboratory Group (Thailand) Co.,Ltd.
104 Phatthanakan 40, Phatthanakan Rd.,
Khwaeng Phatthanakan, Khet Suan Luang,
Bangkok 10250 Thailand

Ambient Temperature : (20 ± 2.5) °C
Relative Humidity : (50 ± 10) %
Barometric Pressure : 760 mmHg
Calibration Procedure : ASTM E 542 - 01

Calibrated by : Natcha Chayyingcheiw

Approved by : 
Approved Signatory

() Unnopphol Harachai
(✓) Srisuda Khamtha
() Sa-ngeunkam Wongsu

Issue Date : 27 February 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Equipment : Burette
Received Date : 23 February 2024
Condition As-Received : New Item
Calibration Date : 27 February 2024
Reference : 2402-0757DSC-1

Cert.No.: 24CG952
Page.: 2 of 2

Condition of this result of calibration

1. Reference Standard Instruments :

Instruments	Model	Serial No.	ID. No.	Certificate No.	Traceability	Due date
1) Balance	XP205DR	1126143764	140RC004	23MM538	TPA	15 Sep 2024
2) Thermo-Hygrograph	THDX-CE	00016540	140EC001	23H1275	TPA	09 June 2024
3) Thermometer	-	0834181	140EC005	23I948	TPA	10 Aug 2024

This certification is traceable to SI Unit

- The certificate is valid only to the item calibrated on date and place of calibration.
- True value is converted to true volume at the standard temperature of 20 °C

Calibration result :

Nominal capacity (mL)	Reading (mL)	Uncertainty ($\pm$ mL)	k Factor
50	50.0032	0.010	2.00

Remark mL = cm³

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k , providing a level of confidence of approximately 95 %.

-o0o-



REVIEW BY Autcharawan S.
APPROVED BY Tanyatorn M.
NEXT CAL. DATE 12 Jan 2025

Certificate of Calibration
ICS-2100: Anion (ID#659)

This certificate is to verify that instrument below are calibrated

by Archemica Lab Co., Ltd.

ICS-2100 S/N: 15010977

AS-HV S/N: 5450A36659

For

ALS Laboratory Group (Thailand) Co., Ltd.



Operator Signature: Nutdanai Date: Jan 12, 2024

(Mr.Nutdanai Laekhwan)

Application Chemist



Metrological Center

SCI ECO Services Company Limited

33/2 Moo 3, T.Banpa, A.Kaengkhoi, Saraburi 18110

Telephone : +66 2 586 5792-4 Fax : +66 2 586 5109

Website : www.scieco.co.th E-Mail : calibrate@scg.co.th

Certificate No. T240977

Page 1 of 5

Certificate of Calibration

Equipment : DIGESTION UNIT

Manufacturer : Environmental Express

Model : TKN100

Serial No. : 2017TKNBC142

Customer Code : BKK_EN0223

ID No. : T6773A4

Customer : ALS Laboratory Group (Thailand) Co.,Ltd.

104 Phatthanakan 40, Phatthanakan Rd., Khwaeng Phatthanakan,

Khet Suan Luang, Bangkok 10250

Customer Location : Wet Chemistry Lab1

Date of Receipt : 15 May 2024

Calibrated By : Sujjar Naknakred (Site Calibration Manager)

Approved By :  Preecha Phisassutthikul (Temperature Calibration Manager)

Date of Issue : 28 MAY 2024

The uncertainties are for a confidence probability of approximately 95%.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standard laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the Metrological Center.



Metrological Center

SCI ECO Services Company Limited

33/2 Moo 3, T.Banpa, A.Kaengkhoi, Saraburi 18110

Telephone : +66 2 586 5792-4 Fax : +66 2 586 5109

Website : www.scieco.co.th E-Mail : calibrate@scg.co.th

Certificate No. T240977

Page 2 of 5

Calibration Report

Equipment : DIGESTION UNIT

Date of Calibration : 24 May 2024

Environment : Temperature : 28.7 - 30.0 °C

Line Voltage : 222.8 - 225.9 V

Relative Humidity : 55 - 65 %RH

Condition of this results of calibration :

1. This equipment was calibrated by insert four standard thermocouples type S into its chamber, the other one thermocouple type T use for ambient temperature measurement. The calibration was done in according to WI-T10.

2. Reference Standard Instrument :

Instrument	Model	Instrument No.	Certificate No.	Due Date
TC	Type S	M20A1-(CH17-CH20)	T240714	23 April 2025
DATA LOGGER	34970A	T149	T240714	23 April 2025

3. This certificate is traceable to :

National Institute of Metrology (Thailand) through Metrological Center (NSC-TISI-TIS 17025 CALIBRATION 0244.)

4. Condition of calibrated item : good

Equipment Description :

Time Constant 2 Hour 2 Minute At 380 °C

Fresh Air Damper ☐ Open ☐ Min ☐ Medium ☐ Max

☐ Close

☒ Not Available

5. Adjustment :

() without adjustment

(X) after adjustment

Approved By : 



Metrological Center

SCI ECO Services Company Limited

33/2 Moo 3, T.Banpa, A.Kaengkhoi, Saraburi 18110


Telephone : +66 2 586 5792-4 Fax : +66 2 586 5109

Website : www.scieco.co.th E-Mail : calibrate@scg.co.th

Certificate No. T240977

Page 3 of 5

Calibration Report

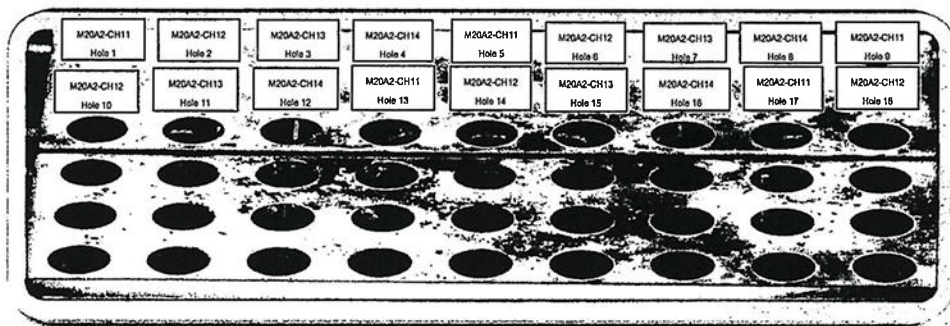

DISPLAY CONTROL (FRONT)

Measurement Results

Cal. Point	Setting	Reading	STD.	Position of Standards at Block								
(°C)	(°C)	(°C)	Reading	M20A2-CH11 Hole 1	M20A2-CH12 Hole 2	M20A2-CH13 Hole 3	M20A2-CH14 Hole 4	M20A2-CH11 Hole 5	M20A2-CH12 Hole 6	M20A2-CH13 Hole 7	M20A2-CH14 Hole 8	M20A2-CH11 Hole 9
380.0	380.0	379.8 - 380.2	Average °C	377.34	380.04	380.62	383.20	380.97	379.00	381.21	378.62	379.36
			Stability ± °C	0.31	0.27	0.31	0.25	0.24	0.37	0.15	0.26	0.23

Cal. Point	Setting	Reading	STD.	Position of Standards at Block								
(°C)	(°C)	(°C)	Reading	M20A2-CH12 Hole 10	M20A2-CH13 Hole 11	M20A2-CH14 Hole 12	M20A2-CH11 Hole 13	M20A2-CH12 Hole 14	M20A2-CH13 Hole 15	M20A2-CH14 Hole 16	M20A2-CH11 Hole 17	M20A2-CH12 Hole 18
380.0	380.0	379.8 - 380.2	Average °C	376.66	381.17	380.28	383.11	383.26	382.04	380.37	381.00	381.61
			Stability ± °C	0.35	0.34	0.35	0.24	0.21	0.23	0.23	0.20	0.22

Approved By.

FM-L13 I08/30-05-57



Metrological Center

SCI ECO Services Company Limited

33/2 Moo 3, T.Banpa, A.Kaengkhoi, Saraburi 18110


Telephone : +66 2 586 5792-4 Fax : +66 2 586 5109

Website : www.scieco.co.th E-Mail : calibrate@scg.co.th

Certificate No. T240977

Page 4 of 5

Calibration Report

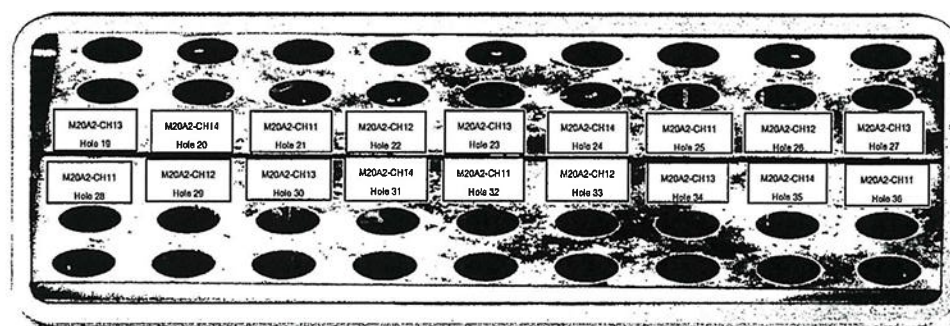

DISPLAY CONTROL (FRONT)

Measurement Results

Cal. Point	Setting	Reading	STD.	Position of Standards at Block								
(°C)	(°C)	(°C)	Reading	M20A2-CH13 Hole 19	M20A2-CH14 Hole 20	M20A2-CH11 Hole 21	M20A2-CH12 Hole 22	M20A2-CH13 Hole 23	M20A2-CH14 Hole 24	M20A2-CH11 Hole 25	M20A2-CH12 Hole 26	M20A2-CH13 Hole 27
380.0	380.0	379.8 - 380.2	Average °C	376.97	383.29	383.26	381.86	380.47	381.36	382.07	379.37	381.87
			Stability ± °C	0.22	0.26	0.22	0.22	0.23	0.22	0.25	0.21	0.22

Cal. Point	Setting	Reading	STD.	Position of Standards at Block								
(°C)	(°C)	(°C)	Reading	M20A2-CH11 Hole 28	M20A2-CH12 Hole 29	M20A2-CH13 Hole 30	M20A2-CH14 Hole 31	M20A2-CH11 Hole 32	M20A2-CH12 Hole 33	M20A2-CH13 Hole 34	M20A2-CH14 Hole 35	M20A2-CH11 Hole 36
380.0	380.0	379.8 - 380.2	Average °C	377.29	381.19	381.69	380.72	381.02	380.42	378.65	380.63	376.92
			Stability ± °C	0.25	0.27	0.27	0.26	0.22	0.28	0.19	0.23	0.22

Approved By.

FM-L13 I08/30-05-57



Metrological Center

SCI ECO Services Company Limited

33/2 Moo 3, T.Banpa, A.Kaengkhoi, Saraburi 18110


Telephone : +66 2 586 5792-4 Fax : +66 2 586 5109

Website : www.scieco.co.th E-Mail : calibrate@scg.co.th

Certificate No. T240977

Page 5 of 5

Calibration Report

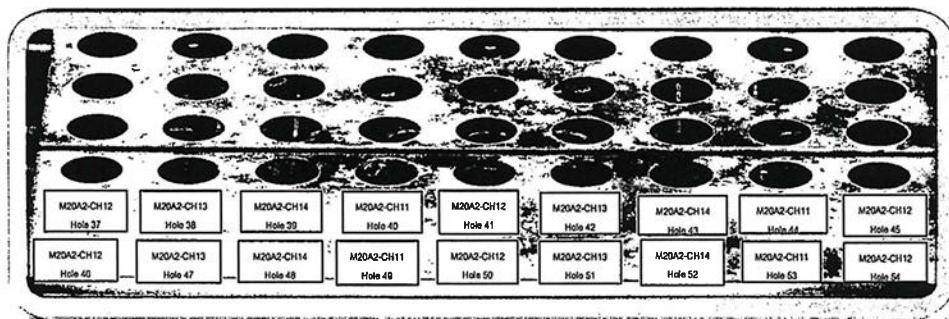

DISPLAY CONTROL (FRONT)

Measurement Results

Cal. Point	Setting	Reading	STD.	Position of Standards at Block								
(°C)	(°C)	(°C)	Reading	M20A2-CH13 Hole 37	M20A2-CH13 Hole 38	M20A2-CH14 Hole 39	M20A2-CH11 Hole 40	M20A2-CH12 Hole 41	M20A2-CH13 Hole 42	M20A2-CH14 Hole 43	M20A2-CH11 Hole 44	M20A2-CH13 Hole 45
380.0	380.0	379.8 - 380.2	Average °C	380.39	378.84	380.87	381.85	379.87	382.97	383.29	379.42	378.50
			Stability ±°C	0.21	0.16	0.19	0.25	0.18	0.19	0.19	0.23	0.42

Cal. Point	Setting	Reading	STD.	Position of Standards at Block								
(°C)	(°C)	(°C)	Reading	M20A2-CH12 Hole 46	M20A2-CH13 Hole 47	M20A2-CH14 Hole 48	M20A2-CH11 Hole 49	M20A2-CH12 Hole 50	M20A2-CH13 Hole 51	M20A2-CH14 Hole 52	M20A2-CH11 Hole 53	M20A2-CH12 Hole 54
380.0	380.0	379.8 - 380.2	Average °C	377.89	381.24	380.93	379.91	380.08	380.61	381.63	383.34	383.04
			Stability ± °C	0.63	0.41	0.24	0.14	0.15	0.15	0.14	0.16	0.38

The expanded uncertainty of temperature measurement was $\pm 1.65^\circ\text{C}$

The calibration result apply only the above calibrated item.

The result of test was found accurate as shown on date and place of test only.

The reported expanded uncertainty is based on a standard uncertainty multiplied by a coverage factor $k=2$, providing a level of confidence of approximately 95 %.

Approved By.

FM-L13 108/30-05-57



บริษัท ดับเบิล เอส ไดแอกโนสติกส์ จำกัด
DOUBLE S DIAGNOSTICS CO., LTD.

4 ซอยอุดมสุข 14 แขวงบางนา เขตบางนา กรุงเทพฯ 10260 โทรศัพท์ (02) 747-7009 โทรสาร (02) 747-7008
4 Soi Udomsuk 14, Bangna, Bangkok 10260 Tel. (02) 747-7009 Fax: (02) 747-7008

Maintenance Plan YEAR : 2024

เดือน	Jan	Feb	Mar	Apr	May	Jun	Jul	Aug	Sep	Oct	Nov	Dec
ดำเนินการ	OK											

Periodical maintenance check list for Konelab

	6M	12M	Notel
1.Diluent-wash tubing change	☒	☐	
2.ISE tubing change	☒	☒	None
3.Syringe check/change		☐	
4.Dispensing check/ change		☐	
5.Waste tubing change when necessary		☐	
6.Lamp check/change	☒	☐	
7.Mixer paddle/paddle change(not Konelab20)		☐	
8.ISE needles check/change		☐	
9.Pump tubing check/ change	☒	☐	
10.Broken/worn out part check /change		☐	
11.Peristaltic pump check /cleaning/ lubrication	☒	☐	
12.Heating check		☐	
13.Cooling check		☐	
14.Dispenser mechanic check/adjustment	☒	☐	
15.Cuvette transfer mechanic check/adjustment	☒	☐	
16.Dispenser movement check/adjustment	☒	☐	
17.Sample/reagent register check/adjustment	☒	☐	
18.Dispensing tubing tightness check	☒	☐	
19.Photometer and optics cleaning/check/adjustment	☒	☐	
20.Workstation PC cleaning if necessary	☒	☐	
21.Mechanic cleaning/lubrication	☒	☐	
22.Instrument cleaning if necessary	☒	☐	
23.Complete analyzer testing with waterblank/QC or sample	☒	☐	
24.Test parameters/Adjustment/config. Save to USB key	☒	☐	
25.UPS Test	☒	☐	

Place: DIS LAB Instrument: K80 Aquakem
Date/Time: 9/1/62 Serial no: 22781
Service done by: 69286 Install date:
Signature of customer: 11/1/62 Date/Time:



บริษัท ดับเบิล เอส ไดแอกโนสติกส์ จำกัด
DOUBLE S DIAGNOSTICS CO., LTD.

4 ซอยอุดมสุข 14 แขวงบางนา เขตบางนา กรุงเทพฯ 10260 โทรศัพท์: (02) 747-7009 โทรสาร: (02) 747-7008
4 Soi Udomsuk 14, Bangna, Bangkok 10260 Tel. (02) 747-7009 Fax: (02) 747-7008

Maintenance Plan YEAR : 2024

เดือน	Jan	Feb	Mar	Apr	May	Jun	Jul	Aug	Sep	Oct	Nov	Dec
Jan												

Periodical maintenance check list for Konelab

	6M	12M	Note
1.Diluent-wash tubing change	☐	☒	
2.ISE tubing change	☒	☒	none
3.Syringe check/change		☒	
4.Dispensing check/ change		☒	
5.Waste tubing change when necessary		☒	
6.Lamp check/change	☐	☒	
7.Mixer paddle/paddle change(not Konelab20)		☒	
8.ISE needles check/change		☒	none
9.Pump tubing check/ chance	☐	☒	
10.Broken/worn out part check /change		☒	
11.Peristaltic pump check /cleaning/ lubrication	☐	☒	
12.Heating check		☒	
13.Cooling check		☒	
14.Dispenser mechanic check/adjustment	☐	☒	
15.Cuvette transfer mechanic check/adjustment	☐	☒	
16.Dispenser movement check/adjustment	☐	☒	
17.Sample/reagent register check/adjustment	☐	☒	
18.Dispensing tubing tightness check	☐	☒	
19.Photometer and optics cleaning/check/adjustment	☐	☒	
20.Workstation PC cleaning if necessary	☐	☒	
21.Mechanic cleaning/lubrication	☐	☒	
22.Instrument cleaning if necessary	☐	☒	
23.Complete analyzer testing with waterblank/QC or sample	☐	☒	
24.Test parameters/Adjustment/config. Save to USB key	☐	☒	
25.UPS Test	☐	☒	

Place: ALS LAB Instrument: K20 Aquaken
Date/Time: 16/8/62 Serial no: 29281
Service done by: B. P. S. Install date:
Signature of customer: อรุณ ธีรวิ Date/Time: 16/08/2024

Accuracy results Aquaken 7.2.AQ2 Page: 1

Laboratory
Analyzer User

8/16/2024 14:53

Performed 8/16/2024
Lot WB34

ACCEPTANCE CRITERIA

	Result	Limit	Warning
Temperature (°C)	37.8	37.0 +/- 1.0	
Dispensing ratio	16.4	14.8 - 17.2	
CV%	0.29	<1.7	
Photometric noise			
Max SD L340 2 (mA)	0.17	<2.0	
Max SD L340 4 (mA)	0.87	<3.0	
Linearity of photometer			
Slope	1.0141	0.94 - 1.06	
Curvature	0.0053	+/- 0.02	
Max bias from linear fit (mA)	4.3	<15.0	
Max delta %	-1.6	+/- 6.0	
Linearity of sample dispensing			
Proport. volume XDISP2 (?l)	2.06	1.96 - 2.16	
Proport. volume XDISP4 (?l)	4.14	3.85 - 4.40	
XDISP2 CV%	1.21	<2.0	
XDISP4 CV%	0.90	<2.0	
XDISP10 CV%	0.68	<2.0	
Needle 0 ?l volume			
Average (A)	0.005	<0.050	
Standard deviation (A)	0.002	<0.005	
Volume (?l)	0.03	<0.32	

OTHER INFORMATION


Dispensing ratio	Photom. noise: SD (mA)
Posit Result (A)	Posit L340 2 L340 4
1 0.1549	1 0.15 0.80
2 0.1549	2 0.17 0.79
3 0.1537	3 0.04 0.65
4 0.1547	4 0.16 0.31
5 0.1547	5 0.11 0.58
6 0.1545	6 0.14 0.87

Linearity of sample dispensing

Test	Absorbance (A)
XDISP2	0.306
XDISP4	0.612
XDISP10	1.471

Linearity of photometer

L340	Target (A)	Meas (A)	Delta (A)	Delta %
1	0.002	0.006	-0.004	-217.7
2	0.486	0.493	-0.007	-1.5
3	1.451	1.469	-0.018	-1.2
4	1.936	1.963	-0.027	-1.4
5	2.415	2.454	-0.039	-1.6

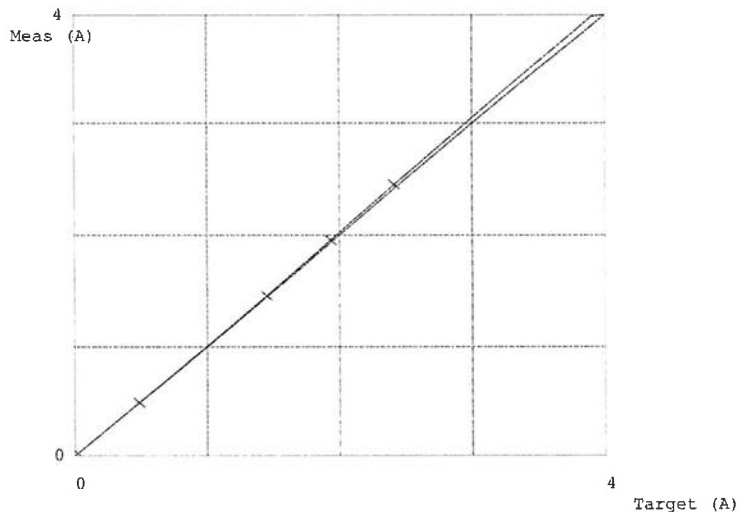

Sartorius (Thailand) Co., Ltd.
129 Rama 9 Road, Huaykwang, Huaykwang, Bangkok 10310
Tel: +66 2643 8361-6, e-mail: service.thailand@sartorius.com



SARTORIUS

NSC-TISI-TIS 17025
CALIBRATION 0426

REVIEW BY Sithichok
APPROVED BY [Signature]
NEXT CAL. DATE 23/01/25

Certificate of Calibration

Model Number : CP3202S Certificate No. : 24BCI0026
Description : Precision Balance Issued Date : Thursday, January 25, 2024
Serial Number : 16704832 Reference No. : 227439
ID No. BKK_ML0001
Manufacturer : Sartorius Page No. : 1 of 2

Customer Name : ALS Laboratory Group (Thailand)Co., Ltd.
104 Phatthanakan 40, Phatthanakan Rd., Khwaeng Phatthanakan, Khut Suan Luang, Bangkok 10250.

Calibrated Place : Balance Room.

Calibrated By : Mr.Chonchai Inthana
Calibration Date : Tuesday, January 23, 2024
Calibration Procedure No. : This calibration was conducted by
Using in-house calibration procedure number (WI-003)
Based on UKAS LAB 14 : 2019

Metrological data : Capacity : 3,200 g Readability : 0.01 g
Reasons for calibration ☐ New Installation ☐ Service / Repaired ☒ Re-calibration/ Maintenance
Ambients Conditions: Temperature : 22.0 °C ± 3.0 °C
Humidity : 60.0 % RH ± 5.0 % RH
Pressure : ±
Equipment Condition: ☒ Good Operate ☐ Fair

Measurement Method UKAS Publication Ref :Lab 14

The measurement uncertainty stated is the expended uncertainty which is obtained from the standard uncertainty multiplied by the coverage factor (k=2) to provide a level of confidence of approximately 95%. It is determined in accordance with the Guide to Expression of Uncertainty in Measurement (GUM). The calibration certificate documents the traceability to National Standards, which realise the unit of measurement according to the International Standard System of Units (SI). Report of Tolerance came form list of Sartorius Metrological Specifications.

Traceability:

Model Number	Description	Traceability	Cetificate No.	Due Date
YCS011-522-00	Sartorius weight set 1mg - 5000g E2,YCS011-522-00	TCS	M2308197S	23-Aug-2025
MHB-382SD	Humidity/Barometer/Temp Lutron MHB-382SD	DKSH	C19231845	23-Aug-2024

This certificate relate and apply this equipment only.

This certificate may not be reproduced other than in full except with the prior written approval of the Verification Operation Division
Sartorius (Thailand) Co., Ltd.

[Signature]

Mr.Chonchai Inthana(Technical Manager)

SOP FM 33 03 February 2022



Sartorius (Thailand) Co., Ltd.

129 Rama 9 Road, Huaykwang, Huaykwang, Bangkok 10310
Tel: +66 2643 8361-6 Fax: +66 2643-8397, e-mail: service.thailand@sartorius.com


SARTORIUS

Certificate of Calibration

Model Number : CP3202S
Description : Precision Balance
Serial Number : 16704832
ID No. : BKK_ML0001
Manufacturer : Sartorius

Certificate No. : 24BCI0026
Issued Date : Thursday, January 25, 2024
Reference No. : 227439
Page No. : 2 of 2

Calibration Results : Without Adjustment

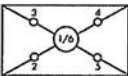Repeatability			Eccentricity (Off-center loading error)																
<i>The reproducibility is the ability of a weighing instrument to display nearly identical readouts under constant test conditions when the same load within a measurement series is placed repeatedly on the weighing pan in the same manner. The standard deviation is used to express reproducibility quantitatively.</i>			<i>The off-center loading error is yielded by the difference between the readout of the load, i.e. 1/3 or 1/4 of maximum capacity, placed in the middle of the weighing pan and between each of four additional measurement points (positions defined according to OIML R76).</i>																
Nominal Value : (Low Load)	200.00	2000.00	Nominal value :	1000	g														
200 g	200.00	2000.00	Tolerance	0.06	g														
Tolerance	200.00	2000.00																	
0.01 g	200.00	2000.00	<table><tr><th colspan="2">Difference</th></tr><tr><td>1</td><td>—</td></tr><tr><td>2</td><td>-0.01</td></tr><tr><td>3</td><td>0.00</td></tr><tr><td>4</td><td>0.00</td></tr><tr><td>5</td><td>0.00</td></tr><tr><td>6</td><td>-</td></tr></table>			Difference		1	—	2	-0.01	3	0.00	4	0.00	5	0.00	6	-
Difference																			
1	—																		
2	-0.01																		
3	0.00																		
4	0.00																		
5	0.00																		
6	-																		
	200.00	1999.99																	
Nominal Value : (High Load)	199.99	2000.00																	
2000 g	200.00	2000.00																	
Tolerance	200.00	2000.00																	
0.01 g	200.00	2000.00																	
	200.00	1999.99																	
Standard Deviation	0.003	0.004																	

Linearity

The linearity, also called linearity error. Describes the deviation of the characteristic curve of a weighing instrument from the linear slope.

Tolerance 0.02 g				
Nominal Value (g)	Conventional Mass Value (g)	Displayed Value (g)	Deviation (g)	Uncertainty (g)
0.5	0.50	0.50	0.00	0.011
1	1.00	1.00	0.00	0.011
10	10.00	10.00	0.00	0.011
100	100.00	100.00	0.00	0.011
200	200.00	200.00	0.00	0.011
500	500.00	500.00	0.00	0.011
1000	1000.00	1000.00	0.00	0.011
1500	1500.00	1500.00	0.00	0.012
2000	2000.00	2000.00	0.00	0.012
3000	3000.00	3000.01	0.01	0.012

End of Report.



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000-29 FAX. 0-2719-9484



Cert. No.: 23TM1408

Page : 1 of 4

Certificate of Calibration

Equipment : Autoclave
Manufacturer : TOMY
Model : SX-700
Serial No. : 48134190
ID No. : BKK_ML0041
Submitted by : ALS Laboratory Group (Thailand) Co., Ltd.
104 Phatthanakan 40, Phatthanakan Rd.,
Khwaeng Phatthanakan, Khet Suan Luang,
Bangkok 10250 Thailand
Location : Media Preparation Room
Received Order : 03 October 2023
Calibration Date : 04 October 2023
Ambient Temperature : (26 ± 10) °C
Relative Humidity : (50 ± 30) %
Calibrated by : Khit Ruttanaprapachai

Approved by :

Approved Signatory

() Pornthippa Tameyakul
(✓) Ponpan Paipim
() Suwit Imjai

Issue Date : 11 October 2023

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Equipment : Autoclave
 Condition As-Received : Used Item
 Reference : 2310-0006OC-6
 Procedure Used :-

Cert. No.: 23TM1408

Page : 2 of 4

Calibration were conducted using in-house calibration procedure CP-OT03 according to direct measurement method with Data Acquisition which connected with Thermocouple Type T

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Data Acquisition	MY57013823	23LM66	TPA	25 Mar 2024

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

4. This result of calibration covers laboratory autoclaves for the sterilization of goods and material which could be infected with organisms categorized as Hazard Group 1, 2 and 3**

(** = Categorization of pathogens according to hazard and categories of containment, second edition, 1990)


It does not cover autoclaves for use with material infect with organisms in Hazard Group 4, for which complete containment and sterilization of infected condensate is considered to be essential.

This result of calibration does not apply to sterilizers or disinfectors used for medical, dental, pharmaceutical or veterinary purposes which are directly concerned with patient care, or those used for fabrics subjected to sterilization which are required to be dry at the end of cycle.

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

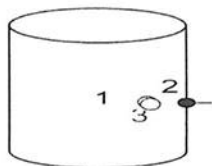

	Environmental		
	(°C)	(%R.H.)	(Volt)
Beginning of Calibration	26	64	221
Finished of Calibration	27	67	222

Position	Description	Ref. Std. ID No.:
1 =	Center of chamber	19-17TC-08
2 =	Temperature sensor	19-17TC-09
3 =	Exhaust port	19-17TC-10

a 1184533



Equipment : Autoclave
 Condition As-Received : Used Item
 Reference : 2310-0006OC-6
 Result of Calibration :- (*) Without Adjustment
 Function of UUC* : Temperature Source

Cert. No.: 23TM1408

Page : 3 of 4

Operating parameter Set : Temperature = 108 °C

Sterilization period = 10 minute

UUC* Setting (°C)	UUC* Reading (°C)	Position	Average* Standard Reading (°C)	Stability (± °C)	Pressure Reading (MPa)	Uncertainty (± °C)	Coverage Factor k
108	108	1	108.352	0.12	0.04	0.90	2
		2	108.263				
		3	108.140				

Operating parameter Set : Temperature = 115 °C

Sterilization period = 20 minute

UUC* Setting (°C)	UUC* Reading (°C)	Position	Average* Standard Reading (°C)	Stability (± °C)	Pressure Reading (MPa)	Uncertainty (± °C)	Coverage Factor k
115	115	1	115.376	0.13	0.08	0.90	2
		2	115.297				
		3	115.157				

Operating parameter Set : Temperature = 118 °C

Sterilization period = 10 minute

UUC* Setting (°C)	UUC* Reading (°C)	Position	Average* Standard Reading (°C)	Stability (± °C)	Pressure Reading (MPa)	Uncertainty (± °C)	Coverage Factor k
118	118	1	118.083	0.11	0.09	0.90	2
		2	118.037				
		3	117.954				

Average* : The average of 30 values in each position.

Stability : One-half of the greatest maximum difference of measured temperature at any one probe.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

a 1184532



Equipment : Autoclave
Condition As-Received : Used Item
Reference : 2310-0006OC-6
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source

Cert. No.: 23TM1408
Page : 4 of 4

Operating parameter Set : Temperature = 121 °C
Sterilization period = 30 minute

UUC* Setting (°C)	UUC* Reading (°C)	Position	Average* Standard Reading (°C)	Stability (± °C)	Pressure Reading (MPa)	Uncertainty (± °C)	Coverage Factor <i>k</i>
121	121	1	121.186	0.17	0.11	0.91	2
		2	121.082				
		3	120.980				

Average* : The average of 30 values in each position.

Stability : One-half of the greatest maximum difference of measured temperature at any one probe.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was excluded stability.

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

a 1184531



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000-29 FAX. 0-2719-9484



Cert. No.: 23TM1366
Page : 1 of 3

Certificate of Calibration

Equipment : Incubator
Manufacturer : Memmert
Model : ICP 750
Serial No. : F816.0004
ID No. : BKK_ML0189
Submitted by : ALS Laboratory Group (Thailand) Co.,Ltd.
104 Phatthanakan 40, Phatthanakan Rd.,
Khwaeng Phatthanakan, Khet Suan Luang ,
Bangkok 10250 Thailand
Location : Incubation & Microbiological Reading
Received Order : 03 October 2023
Calibration Date : 03 October 2023
Ambient Temperature : (26 ± 10) °C
Relative Humidity : (50 ± 30) %
Calibrated by : Man Pattanapongpaiboon

Approved by :
() Pornthippa Tameyakul
(✓) Ponpan Paipim
() Suwit Imjai

Approved Signatory

Issue Date : 11 October 2023

The Uncertainties are for a confidence probability of approximately 95 %

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

A 0053222



Equipment : Incubator
Condition As-Received : Used Item
Reference : 2310-0006OC-3

Cert. No.: 23TM1366

Page : 2 of 3

Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD).

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-


- | Instrument | Serial No. | Cert. No. | Traceable | Due Date |
|----------------------|------------|-----------|-----------|-------------|
| 1) Data Acquisition | MY49023932 | 23LM122 | TPA | 26 Jul 2024 |
2. This certificate is valid only to the item calibrated on date and place of calibration.
3. This certification is traceable to the International System of Unit.

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Close

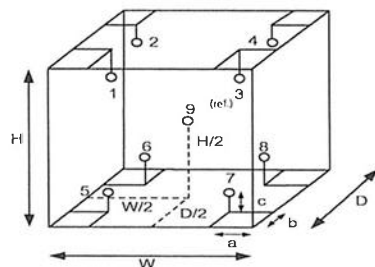

Probe Installation Details :

a = 10 cm
b = 10 cm
c = 10 cm

Dimension of Chamber :

D = 0.50 m
W = 1.1 m
H = 0.80 m
Capacity = 0.44 m³

Environment during calibration		
	Beginning	Finished
Temp. (°C)	24	24
REL.Humid. (%)	54	56
AC Supply (Volt)	222	224

Position :	Ref. Std. ID No.:
1	19-16RTD-01
2	19-16RTD-02
3	19-16RTD-03
4	19-16RTD-04
5	19-16RTD-05
6	19-16RTD-06
7	21-16RTD-07
8	19-16RTD-08
9 (ref.)	19-16RTD-09



Equipment : Incubator
Condition As-Received : Used Item
Reference : 2310-0006OC-3

Cert. No.: 23TM1366

Page : 3 of 3

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Close

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor k
37.0	37.0	37.0	0.026	0.13	0.18	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	
37.0	36.994	36.907	37.045	37.002	37.007	37.002	37.008	36.995	36.923	0.30

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert. No.: 24TM1335
Page : 1 of 3

Equipment : Incubator
Manufacturer : Memmert
Model : ICP 750
Serial No. : F816.0004
ID No. : BKK_ML0189

Submitted by : ALS Laboratory Group (Thailand) Co.,Ltd.
104 Phatthanakan 40, Phatthanakan Rd.,
Khwaeng Phatthanakan, Khet Suan Luang ,
Bangkok 10250 Thailand

Location : Incubation & Microbiological Reading

Received Order : 26 September 2024
Calibration Date : 26 September 2024
Ambient Temperature : (26 ± 10) °C
Relative Humidity : (50 ± 30) %
AC Line Voltage : (220 ± 22) V

Calibrated by : Preecha Hlahib

Approved by :

() Ponpan Paipim
() Suwit Imjai
(✓) Kunchit Promprat

Issue Date : 09 October 2024

REVIEW BY Sithichok T.

APPROVED BY [Signature]

NEXT CAL DATE 26/03/26



Equipment : Incubator
Condition As-Received : Used Item
Reference : 2409-0695OC-2
Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD).
The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Data Acquisition	MY57013711	24LM115	TPA	13 Jul 2025

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.


Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Not Available

Environment during calibration		
	Beginning	Finished
Temp. (°C)	22	22
REL.Humid. (%)	65	59
AC Supply (Volt)	220	221

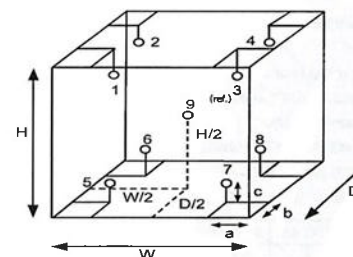

Probe Installation Details :

a = 5.0 cm
b = 5.0 cm
c = 5.0 cm

Dimension of Chamber :

D = 0.60 m
W = 1.0 m
H = 1.2 m
Capacity = 0.75 m³

Position :	Ref. Std. ID No.:
1	18-18RTD-01
2	24-18RTD-02
3	18-18RTD-03
4	18-18RTD-04
5	18-18RTD-05
6	23-18RTD-06
7	18-18RTD-07
8	22-18RTD-08
9 (ref.)	24-18RTD-09

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Equipment : Incubator
 Condition As-Received : Used Item
 Reference : 2409-0695OC-2
 Result of Calibration :- (*) Without Adjustment
 Function of UUC* : Temperature Source
 Fresh air setting : Not Available

Cert. No.: 24TM1335
 Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor <i>k</i>
37.0	37.0	37.0	0.12	0.17	0.30	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	
37.0	36.898	36.906	36.974	36.911	36.943	36.936	36.978	36.958	36.850	0.30

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
 CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
 534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
 TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert. No.: 24TM667
 Page : 1 of 3


Equipment : Hot Air Oven
 Manufacturer : Binder
 Model : ED 240/E2
 Serial No. : 00-15533
 ID No. : BKK_ML0013

Submitted by : ALS Laboratory Group (Thailand) Co.,Ltd.
 104 Phatthanakan 40, Phatthanakan Rd.,
 Khwaeng Phatthanakan, Khet Suan Luang,
 Bangkok 10250 Thailand

Location : Media Preparation Room

Received Order : 23 April 2024
 Calibration Date : 23 April 2024
 Ambient Temperature : (26 ± 10) °C
 Relative Humidity : (50 ± 30) %

Calibrated by : Tawatchai Pama

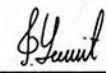Approved by : 
 Approved Signatory

() Ponpan Paipim
 (✓) Suwit Imjai
 () Kunchit Promprat

Issue Date : 26 April 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Equipment : Hot Air Oven
Condition As-Received : Used Item
Reference : 2404-0439OC-8

Cert. No.: 24TM667
Page : 2 of 3

Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Thermocouple Type T.

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Data Acquisition	MY49001451	24LM44	TPA	17 Mar 2025

2. This certificate is valid only to the item calibrated on date and place of calibration.


3. This certification is traceable to the International System of Unit.

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Close

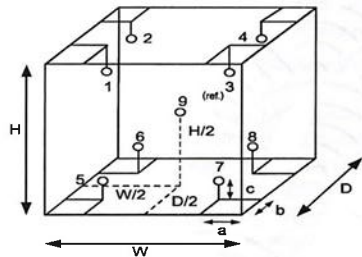

Environment during calibration		
	Beginning	Finished
Temp. (°C)	24	23
REL.Humid. (%)	65	65
AC Supply (Volt)	223	222

Position :	Ref. Std. ID No.:
1	24-19TC-01
2	24-19TC-02
3	24-19TC-03
4	24-19TC-04
5	24-19TC-05
6	24-19TC-06
7	24-19TC-07
8	24-19TC-08
9 (ref.)	24-19TC-09

Probe Installation Details :

a = 10 cm
b = 10 cm
c = 10 cm

Dimension of Chamber :

D = 0.50 m
W = 0.80 m
H = 0.60 m
Capacity = 0.24 m³



Equipment : Hot Air Oven
Condition As-Received : Used Item
Reference : 2404-0439OC-8
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Close

Cert. No.: 24TM667
Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor k
180	180	180	0.64	2.7	3.7	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	
180	181.009	181.511	180.922	181.359	181.217	183.659	181.664	181.986	181.474	1.5

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

ภาคผนวก จ

สำเนาหนังสือใบอนุญาตขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน



ที่อก ๐๓๑๐(๑)/ ๑ ๖ ๑ ๖ ๘

กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๒ ๐ พฤศจิกายน ๒๕๖๖

เรื่อง ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

เรียน กรรมการผู้จัดการ บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๔ สิงหาคม ๒๕๖๖

- สิ่งที่ส่งมาด้วย ๑. รายชื่อผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๑ แผ่น
๒. รายชื่อเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๕ แผ่น
๓. ขอบข่ายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๓๑ แผ่น

ตามหนังสือที่อ้างถึง บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด ขอต่ออายุหนังสือ
รับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ๖-๒๐๔ สถานที่ตั้งเลขที่ ๑๐๔ ซอยพัฒนาการ ๔๐
ถนนพัฒนาการ แขวงพัฒนาการ เขตสวนหลวง กรุงเทพมหานคร ต่อกรมโรงงานอุตสาหกรรม นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด
ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน โดยมีองค์ประกอบดังนี้

- ก. ผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๖ ราย ตามสิ่งที่ส่งมาด้วย ๑
ข. เจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ ๑๘๑ ราย ตามสิ่งที่ส่งมาด้วย ๒
ค. ขอบข่ายสารมลพิษที่ได้รับขึ้นทะเบียนให้วิเคราะห์ในน้ำเสีย น้ำใต้ดิน อากาศเสีย สิ่งปฏิกูล
หรือวัสดุที่ไม่ใช่แล้ว และดิน ตามสิ่งที่ส่งมาด้วย ๓

หนังสือฉบับนี้จะหมดอายุในวันที่ ๒ กันยายน ๒๕๖๙ หากประสงค์จะต่ออายุหนังสือ
รับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ให้ยื่นคำขอต่ออายุพร้อมเอกสารประกอบคำขอต่อ
กรมโรงงานอุตสาหกรรม ภายใน ๓๐ วัน ก่อนวันสิ้นอายุของหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
ทั้งนี้ สามารถยื่นคำขอผ่านระบบอิเล็กทรอนิกส์ได้ที่หน้าเว็บไซต์กรมโรงงานอุตสาหกรรม

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

(นายจิระ จันทรเจิด)

นักวิทยาศาสตร์เชี่ยวชาญ วิชาการจัดการแผน
ผู้อำนวยการกองวิจัยและเตือนภัยมลพิษโรงงาน
ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม

กองวิจัยและเตือนภัยมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕

โทรสาร ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๔๔

ไปรษณีย์อิเล็กทรอนิกส์ saraban@diw.mail.go.th



"อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว"



เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด เลขทะเบียน ๖-๒๐๔

ที่อก ๐๓๑๐(๑)/ ๑ ๖ ๑ ๖ ๘ ลงวันที่ ๒ ๐ พฤศจิกายน ๒๕๖๖

ก. ผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๖ ราย

- | | |
|--------------------------------|----------------------------|
| ๑) นางสาวยุพาพร จันทร์เปล่ง | ทะเบียนเลขที่ ๖-๒๐๔-ค-๐๐๐๑ |
| ๒) นางสาวชัชชัย โกมารกุล ณ นคร | ทะเบียนเลขที่ ๖-๒๐๔-ค-๐๐๐๒ |
| ๓) นายศรายุทธ จิตรานนท์ | ทะเบียนเลขที่ ๖-๒๐๔-ค-๐๐๐๓ |
| ๔) นางสาวกนกกร เอนก | ทะเบียนเลขที่ ๖-๒๐๔-ค-๐๐๐๔ |
| ๕) นายสุริยา สอนแก้ว | ทะเบียนเลขที่ ๖-๒๐๔-ค-๐๐๐๕ |
| ๖) นายวิชาญ ชุมพรี | ทะเบียนเลขที่ ๖-๒๐๔-ค-๐๐๐๖ |

เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด เลขทะเบียน ว-๒๐๔

ที่ ออก ๐๓๑๐(๑)/ ๑๖๑๖๘ ลงวันที่ ๒๐ พฤศจิกายน ๒๕๖๖

ข. เจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๑๘๑ ราย

๑) นายกาจบัณฑิต กิตติศุภวณิชย์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๐๑
๒) นายภัทรพล สว่างใจธรรม	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๐๒
๓) นายนราธิป เทือกชัยคำ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๐๓
๔) นายศิริโชค พงษ์ประสม	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๐๔
๕) นายณัฐวุฒิ ดั่งแสง	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๐๕
๖) นางสาวจินดา โชกุลธรรม	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๐๖
๗) นางสาวสวรินทร์ น้อยเสงี่ยม	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๐๗
๘) นางสาวชนัญญาญจน์ อัมขม	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๐๘
๙) นางสาวนรินทร์ สายแสง	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๐๙
๑๐) นางสาวนันทิต สมบูรณ์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๑๐
๑๑) นางสาวศรียา เฉลิมธำรงค์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๑๑
๑๒) นางสาวอัญญธร มงคลจิรวุฒิ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๑๒
๑๓) นางสาวศิริลักษณ์ บุญนาค	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๑๓
๑๔) นายนพพงศ์ จันทร์พันธุ์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๑๔
๑๕) นายนรเศรษฐ์ โกมลย์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๑๕
๑๖) นายธันวา จรรย์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๑๖
๑๗) นางสาวเกศรินทร์ แก้วมัน	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๑๗
๑๘) นางสาวสุวิมล ชัยเรืองวุฒิ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๑๘
๑๙) นางสาวสุชาดา ธรรมถาวร	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๑๙
๒๐) นางสาวเปมิกา ชัยเดชธนกุล	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๒๐
๒๑) นางสาวศศิธร หมูสวัสดิ์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๒๑
๒๒) นางสาวเสาวลักษณ์ ภูนาอำพร	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๒๒
๒๓) นายอภิสิทธิ์ สิงหา	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๒๓
๒๔) นายศักดิ์สิทธิ์ ไพศาลพิสุทธิ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๒๔
๒๕) ว่าที่ร้อยตรีหญิง พรรณิกา ขำเจริญ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๒๕
๒๖) นางจิตตา คำแก้ว	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๒๖
๒๗) นางสาวอรรพรรณ รักยง	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๒๗
๒๘) นางสาวนพรัตน์ แยมกรานต์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๒๘
๒๙) นายจุลเดช วารินทร์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๒๙
๓๐) นางสาวดาญรัตน์ ร้องคำ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๓๐
๓๑) นายพรมมี ศรีปัตเนตร	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๓๑
๓๒) นายอุทิศ อุณสิม	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๓๒
๓๓) ว่าที่ร้อยตรี เฉลิมเกียรติ อมรศรีเสริม	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๓๓
๓๔) นางสาววริยา สร้างนา	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๓๔
๓๕) นายอนุพงษ์ รัตนศรีประเสริฐ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๓๕

วิมล

๓๖) นางสาวจุฑารัตน์...

๓๖) นางสาวจุฑารัตน์ โอนสันเทียะ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๓๖
๓๗) นางสาวจาวรรณ พิมพ์อภิกฤติยา	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๓๗
๓๘) นางสาวปรารถน์ทิพย์ กิจไพศาลศักดิ์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๓๘
๓๙) นางสาวเตือนใจ ทางกลาง	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๓๙
๔๐) นางสาวจิราพร ศิริเวช	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๔๐
๔๑) นายวรกร ผูกักษ์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๔๑
๔๒) นายทง วิริยะสทกิจ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๔๒
๔๓) นายธนิศ เจนจบ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๔๓
๔๔) นายคณิศร ขำเพชร	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๔๔
๔๕) นายภูวิช พรหมสะอาด	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๔๕
๔๖) นายธนเดช โกคาพิพัฒน์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๔๖
๔๗) นายชวฤทธิ์ วงษ์จันทร์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๔๗
๔๘) นายอาทิตย์ ศรีเสน	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๔๘
๔๙) นายเจษฎินทร คงศักดิ์ไทย	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๔๙
๕๐) นายจรัส บุญยั้ง	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๕๐
๕๑) นายธนาณัติ เอนก	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๕๑
๕๒) นายอภิวัฒน์ ทุมพู	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๕๒
๕๓) นางสาวสุภาวัญ มาก	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๕๓
๕๔) นางสาวทัศนพร ขวาลสมบูรณ์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๕๔
๕๕) นางสาวอติมา บุญเพ็ง	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๕๕
๕๖) นางสาวภาณุมาศ นามวัฒน์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๕๖
๕๗) นางสาวอุไรรัตน์ ทิงสร้างแป้น	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๕๗
๕๘) นายธีรวัฒน์ ปางสุข	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๕๘
๕๙) นายอิทธิพล ยะโส	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๕๙
๖๐) นายประพจน์ วรรณชูชัย	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๐
๖๑) นายชยธร พงทิพย์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๑
๖๒) นางสาวกนกวรรณ จันทบาล	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๒
๖๓) นายสิทธิโชค ธงเงิน	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๓
๖๔) นางศิวารณ ใจบุญ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๔
๖๕) นางสาวพรรณธิดา พุ่มคง	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๕
๖๖) นายนวภัทร ศรีวิริยะ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๖
๖๗) นายสุวิชา ทองอ่อน	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๗
๖๘) นายวิญญู บุญตะนัย	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๘
๖๙) นายสมบูรณ์ บุตรจันทร์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๙
๗๐) นายวิรัตน์ ไชยชนะรา	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๗๐
๗๑) นายนฤพนธ์ เพิ่มพูน	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๗๑
๗๒) นายจิรณัฐ ขาวละออ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๗๒
๗๓) นายอัสรี นามบุรี	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๗๓
๗๔) นายอัครเวศ จ่อสา	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๗๔

วิมล

๗๕) นายประเสริฐ...

๗๕) นายประเสริฐ สุระขันธ
๗๖) นายบุญล จันทรเนียม
๗๗) นายพรพงษ์ ทองคุณปรีดา
๗๘) นายณฤพล ทองนุช
๗๙) นายอนุวัฒน์ ม่วงแพ
๘๐) นายเจตศราวุฒิ ปิตตะมะ
๘๑) นายกฤษณะ สายวรรณ
๘๒) นายพิชัย บุญยงค์
๘๓) นายภาณุพงศ์ โยมวงศ์
๘๔) นายสามารถ คุ้มปลี
๘๕) นายสัณชัย โกศรินาม
๘๖) นายณัฐวุฒิ ศรีประเสริฐ
๘๗) นายชวลิตชัย นาคพนม
๘๘) นายพศธร ชัยทิพย์
๘๙) นายสิทธิโชค ทาสีดา
๙๐) นายธนากร อินสุตา
๙๑) นางสาววรรณิษา ขาติวันชัย
๙๒) นางสาวพิมพ์ตะวัน มินากุล
๙๓) นางสาวเพชรรัตน์ สิงห์สมบุญ
๙๔) นางสาวชญาณีน พรหมจันทร์
๙๕) นายกิริติ ทวีราช
๙๖) นายจักริน หมั่นวิชา
๙๗) นายฉัตรชัย สุขเปี้ย
๙๘) นายณรรนท ต๊ะทองคำ
๙๙) นายดุลยพล สนนอก
๑๐๐) นายทักษ์ดนัย อุบลศรี
๑๐๑) นายธนศร นามะกฤษณา
๑๐๒) นายธิตีพงศ์ บัวแดง
๑๐๓) นายณนทชัย อุปถัมภ์
๑๐๔) นายณัฐพล คุณสุทธิ
๑๐๕) นายณัฏฐวัฒน์ สาริน
๑๐๖) นายปิยะนัฐ พลมะศรี
๑๐๗) นายพงศ์สิริ โสมเขียว
๑๐๘) นายพรพัฒน์ กำคำ
๑๐๙) นายภาณุพงศ์ มานิตย์
๑๑๐) นายมงคล ผลาทิพย์
๑๑๑) นายสิรินท ทองอัน
๑๑๒) นายอนชา ทันสมัย
๑๑๓) นายอดิศักดิ์ ผมไผ

ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๗๕
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๗๖
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๗๗
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๗๘
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๗๙
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๘๐
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๘๑
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๘๒
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๘๓
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๘๔
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๘๕
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๘๖
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๘๗
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๘๘
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๘๙
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๙๐
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๙๑
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๙๒
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๙๓
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๙๔
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๙๕
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๙๖
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๙๗
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๙๘
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๐๙๙
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๐๐
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๐๑
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๐๒
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๐๓
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๐๔
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๐๕
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๐๖
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๐๗
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๐๘
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๐๙
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๑๐
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๑๑
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๑๒
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๑๓

วิภา

๑๑๔) นายอนันตชัย...

๑๑๔) นายอนันตชัย วิสม
๑๑๕) นายวรุต ดีนัก
๑๑๖) นายแสงตะวัน นทะสีต
๑๑๗) นายยุทธพงศ์ รัตนะ
๑๑๘) นายชัยณัฐมิ ไชยชนะ
๑๑๙) นายวิศรุต ศรีธรรมมา
๑๒๐) นายณนทกร เผือกผ่อง
๑๒๑) นายกำชัย สุทธะ
๑๒๒) นางสาวณัฐภรณ์ บุญตะนัย
๑๒๓) นางสาวพัชรินทร์ แสนสร้อย
๑๒๔) นายไพโรจน์ เปี่ยมพิมาย
๑๒๕) นางสาวศุภมาศ ทองมาก
๑๒๖) นางสาวลลิตา จิตรสว่าง
๑๒๗) นางสาวไมพร เลิกภูเขียว
๑๒๘) นางสาวกฤติมาพร คำมีแก่น
๑๒๙) นางสาวสกุศลรัตน์ ภาควง
๑๓๐) นางสาวไพรินทร์ ศรีรูป
๑๓๑) นางสาวทิพนตร พุฒปัญญา
๑๓๒) นางสาวสาธิตา ปานทอง
๑๓๓) นางสาวอริสา ทองนวล
๑๓๔) นางสาวอรยา คำคล้อง
๑๓๕) นางสาวบุษมากร สุนทรสนาน
๑๓๖) นางสาวอัญชลี คำจันทร์
๑๓๗) นายบุญฤทธิ์ เอี่ยมเทศ
๑๓๘) นางสาวศุภรดา ปันมยุรา
๑๓๙) นางสาวพาดิ คุณนาน
๑๔๐) นางสาวจิราเจต พองดา
๑๔๑) นางสาวอารยา มีชัย
๑๔๒) นางสาววิษุตา นาคผจญ
๑๔๓) นางสาวนันทิยา จันทะสุน
๑๔๔) นายกิตติพงศ์ แซ่ลี
๑๔๕) นายอนุวัติ ภูถวิล
๑๔๖) นายธีรพล แสงทอง
๑๔๗) นายศักดิ์พัฒน์ บุญมัน
๑๔๘) นายจิตติวัตร เอมอุไร
๑๔๙) นายชัยณรงค์ ศรีบุรินทร์
๑๕๐) นางสาวอัจฉราวรรณ สวนสนอง
๑๕๑) นางสาวณัฐพร สิงหา
๑๕๒) นายกัมเรศ แหยมโต

ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๑๔
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๑๕
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๑๖
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๑๗
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๑๘
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๑๙
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๒๐
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๒๑
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๒๒
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๒๓
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๒๔
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๒๕
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๒๖
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๒๗
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๒๘
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๒๙
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๓๐
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๓๑
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๓๒
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๓๓
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๓๔
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๓๕
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๓๖
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๓๗
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๓๘
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๓๙
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๔๐
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๔๑
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๔๒
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๔๓
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๔๔
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๔๕
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๔๖
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๔๗
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๔๘
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๔๙
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๕๐
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๕๑
ทะเบียนเลขที่ ๖-๒๐๔-จ-๐๑๕๒

วิภา

๑๕๓) นางสาวอุบล...

- ๕ -

๑๕๓) นางสาวอุบล เคิกศิริ
๑๕๔) นางสาวมนรัตน์ ทองบุตร
๑๕๕) นายภาคภูมิ แทนไทย
๑๕๖) นางสาวสุภาณัฐ เมลล์พวง
๑๕๗) นางสาวพรทิพา สาตาขนม
๑๕๘) นายเอกวิทย์ วันทะนา
๑๕๙) นายไตรมณฑล ทิพย์วรรณ
๑๖๐) นายจิรเมธ ประเสริฐศิริพงษ์
๑๖๑) นายจิรายุส เกษมสุข
๑๖๒) นายจิรศักดิ์ ศรีวิชัย
๑๖๓) นายณัฐฤกษ์ สหพานแก้ว
๑๖๔) นายบุญธรรมศักดิ์ ปะที
๑๖๕) นายปิ่นณวิญญ์ เสมอทรัพย์
๑๖๖) นายพิษณุพงษ์ ไชยา
๑๖๗) นายภัทรพงษ์ มณฑาทอง
๑๖๘) นายสันต์ ตรีนกุล
๑๖๙) นายภาณุเดช เพชรอุต
๑๗๐) นายอนุกุล วิเศษแสง
๑๗๑) นายภัทรพงษ์ มีสุข
๑๗๒) นางสาวนุชรี สิละทีป
๑๗๓) นางสาวสุภาวดี โกศรินาม
๑๗๔) นางสาวอรณิษ เทียนดำ
๑๗๕) นางสาวพรเพ็ญ ขอบสอน
๑๗๖) นางสาววันวิสา ขอนพิกุล
๑๗๗) นางสาวอรรพรรณ เถาว์ทอง
๑๗๘) นางสาวอัยยิดิน เมอร์วินณ์
๑๗๙) นางสาววิสา คุ้มครอง
๑๘๐) นายวุฒิกร ศิริวรรณ
๑๘๑) นางสาวจางวรรณ กระจำพันธุ์

ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๕๓
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๕๔
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๕๕
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๕๖
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๕๗
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๕๘
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๕๙
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๐
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๑
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๒
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๓
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๔
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๕
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๖
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๗
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๘
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๙
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๗๐
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๗๑
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๗๒
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๗๓
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๗๔
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๗๕
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๗๖
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๗๗
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๗๘
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๗๙
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๐
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๑

วิมล

เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด เลขทะเบียน ว-๒๐๔
ที่ อก ๐๓๑๐(๑)/ ๑ ๖ ๑ ๖ ๘ ลงวันที่ ๒๐ พฤศจิกายน ๒๕๖๖

ค. ขอบข่ายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๓๗๔ รายการ

น้ำเสีย จำนวน 60 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Aldicarb	High-Performance Liquid Chromatographic Method ⁽⁴⁾
2	Aldicarb Sulfone	High-Performance Liquid Chromatographic Method ⁽⁴⁾
3	Aldicarb Sulfoxide	High-Performance Liquid Chromatographic Method ⁽⁴⁾
4	Aldrin	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
5	Arsenic	1) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ⁽⁴⁾
6	Barium	1) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ⁽⁴⁾
7	α-BHC	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
8	β-BHC	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
9	δ-BHC	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
10	γ-BHC	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
11	Biochemical Oxygen Demand	1) 5-Day BOD Test, Azide Modification Method ⁽⁴⁾ 2) 5-Day BOD Test, Membrane Electrode Method ⁽⁴⁾
12	Carbaryl	High-Performance Liquid Chromatographic Method ⁽⁴⁾
13	Carbofuran	High-Performance Liquid Chromatographic Method ⁽⁴⁾
14	Cadmium	1) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ⁽⁴⁾
15	Chemical Oxygen Demand	1) Closed Reflux, Colorimetric Method ⁽⁴⁾ 2) Closed Reflux, Titrimetric Method ⁽⁴⁾
16	Chlordane	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
17	Chromium	1) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ⁽⁴⁾
18	Color	ADMI Weighted-Ordinate Spectrophotometric Method ⁽⁴⁾

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
19	Copper	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
20	Cyanide	Distillation, Colorimetric Method ^[4]
21	2,4'-DDD	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
22	4,4'-DDD	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
23	2,4'-DDE	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
24	4,4'-DDE	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
25	2,4'-DDT	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
26	4,4'-DDT	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
27	Dieldrin	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
28	Endosulfan Sulfate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
29	Endosulfan I	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
30	Endosulfan II	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
31	Endrin	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
32	Endrin Aldehyde	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
33	Formaldehyde	Distillation, Colorimetric Method ^[3]
34	Free Chlorine	1) DPD Ferrous Titrimetric Method ^[4] 2) DPD Colorimetric Method ^[4]
35	Heptachlor	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
36	Heptachlor Epoxide	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
37	Hexavalent Chromium	Colorimetric Method ^[4]
38	3-Hydroxycarbofuran	High-Performance Liquid Chromatographic Method ^[4]
39	Lead	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
40	Manganese	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
41	Mercury	1) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass spectrometric Method ^[4]
42	Methiocarb	High-Performance Liquid Chromatographic Method ^[4]
43	Methoxychlor	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
44	Methomyl	High-Performance Liquid Chromatographic Method ^[4]
45	Nickel	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
46	Oil & Grease	1) Liquid-Liquid, Partition-Gravimetric Method ^[4] 2) Soxhlet Extraction Method ^[4]
47	Oxamyl	High-Performance Liquid Chromatographic Method ^[4]
48	Propoxur	High-Performance Liquid Chromatographic Method ^[4]
49	pH	Electrometric Method ^[4]
50	Phenols	1) Distillation, Chloroform Extraction Method ^[4] 2) Distillation, Direct Photometric Method ^[4]
51	Selenium	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
52	Sulfide	Iodometric Method ^[4]
53	Temperature	Laboratory and Field Methods ^[4]
54	Total Dissolved Solids	Dried at 180 °C ^[4]
55	Total Kjeldahl Nitrogen	Semi-Micro Kjeldahl Method ^[4]
56	Total Phosphorous	Digestion, Colorimetric Method ^[4]
57	Total Suspended Solids	Dried from 103-105 °C ^[4]
58	Toxaphene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
59	Trivalent Chromium	1) Digestion, Inductively Coupled Plasma Method; Colorimetric Method; Calculation ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method; Colorimetric Method; Calculation ^[4]
60	Zinc	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[4]

น้ำใต้ดิน จำนวน 126 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Acenaphthene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
2	Acetone	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
3	Aldrin	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
4	Anthracene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
5	Antimony	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
6	Arsenic	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
7	Atrazine	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
8	Barium	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
9	Benz(a)anthracene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
10	Benzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
11	Benzo(b)fluoranthene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
12	Benzo(k)fluoranthene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
13	Benzoic Acid	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
14	Benzo(a)pyrene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
15	Benzo(g,h,i)perylene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
16	Beryllium	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
17	Bis(2-chloroethyl)ether	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]

18 Bis(2-ethylhexyl)phthalate...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
18	Bis(2-ethylhexyl)phthalate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
19	Bromodichloromethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
20	Bromoform	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
21	Butanol	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
22	Butyl benzyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
23	Cadmium	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
24	Carbazole	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
25	Carbon disulfide	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
26	Carbon tetrachloride	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
27	Chlordane	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
28	p-Chloroaniline	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
29	Chlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
30	Chlorodibromomethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
31	Chloroform	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
32	2-Chlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
33	Chromium	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
34	Chromium (III)	1) Digestion, Inductively Coupled Plasma Method; Colorimetric Method; Calculation ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method; Colorimetric Method; Calculation ^[4]
35	Chromium (VI)	Colorimetric Method ^[4]

36 Chrysene...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
36	Chrysene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
37	Cyanide	Distillation, Colorimetric Method ⁽⁴⁾
38	2,4-D	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
39	DDD	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
40	DDE	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
41	DDT	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
42	Dibenz(a,h)anthracene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
43	Di-n-Butyl Phthalate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
44	1,2-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
45	1,3-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
46	1,4-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
47	3,3-Dichlorobenzidine	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
48	1,1-Dichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
49	1,2-Dichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
50	1,1-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
51	cis-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
52	trans-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
53	2,4-Dichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
54	1,2-Dichloropropane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
55	1,3-Dichloropropane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
56	1,3-Dichloropropene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
57	Dieldrin	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
58	Diethyl Phthalate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
59	2,4-Dimethylphenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
60	2,4-Dinitrophenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
61	2,4-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
62	2,6-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
63	Di-n-octyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
64	Endosulfan	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
65	Endrin	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
66	Ethylbenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
67	Fluoranthene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
68	Fluorene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
69	Heptachlor	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
70	Heptachlor epoxide	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
71	Hexachlorobenzene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
72	Hexachloro-1,3-butadiene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
73	n-Hexane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
74	α-HCH	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾
75	β-HCH	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ⁽⁴⁾

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
76	γ -HCH	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
77	Hexachlorocyclopentadiene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
78	Hexachloroethane	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
79	Indeno(1,2,3-cd)pyrene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
80	Isophorone	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
81	Lead	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
82	Manganese	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
83	Mercury	1) Digestion, Cold Vapor Atomic Absorption Spectrometric Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
84	Methanol	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
85	Methoxychlor	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
86	Methyl bromide	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
87	Methylene chloride	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
88	2-Methylphenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
89	2-Methylnaphthalene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
90	Methyl tert-butyl Ether	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
91	Naphthalene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
92	Nickel	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
93	Nitrobenzene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
94	N-Nitrosodiphenylamine	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
95	N-Nitrosodi-n-Propylamine	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
96	Polychlorinated Biphenyls - PCB 1016 - PCB 1221 - PCB 1232 - PCB 1242 - PCB 1248 - PCB 1254 - PCB 1260	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
97	Pentachlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
98	pH	Electrometric Method ^[4]
99	Phenanthrene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
100	Phenol	1) Distillation, Chloroform Extraction Method ^[4] 2) Distillation, Direct Photometric Method ^[4] 3) Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
101	Pyrene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
102	Selenium	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
103	Silver	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
104	Styrene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
105	1,1,2,2-Tetrachloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
106	Tetrachloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
107	Toluene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
108	Toxaphene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
109	TPH (C ₅ -C ₈)	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[14,25]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
110	TPH (C ₈ -C ₁₆)	Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[9,22]
111	TPH (C ₁₆ -C ₃₅)	Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[9,22]
112	1,2,4-Trichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
113	1,1,1-Trichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
114	1,1,2-Trichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
115	Trichloroethylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
116	2,4,5-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
117	2,4,6-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
118	1,3,5-Trimethylbenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
119	Vanadium	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[4]
120	Vinyl acetate	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
121	Vinyl chloride	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
122	m-Xylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
123	o-Xylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
124	p-Xylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
125	Xylene (Total)	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
126	Zinc	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[4]

อากาศเสีย (ปล่องระบาย) จำนวน 28 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Antimony	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
2	Arsenic	2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5] 1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
3	Beryllium	2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5] 1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
4	Cadmium	2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5] 1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
5	Carbon Monoxide	1) Instrumental Analyzer Method ^[5]
6	Chlorine	2) Sampling Bag Non-Dispersive Infrared Method ^[5] 1) Absorption Sampling, Ion Chromatographic Method ^[5]
7	Chromium	2) Isokinetic Sampling, Ion Chromatographic Method ^[5] 1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
8	Cobalt	2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5] 1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
9	Copper	2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5] 1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
10	Cresol	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
11	Dioxins	Adsorption Sampling, Gas Chromatographic Method ^[5]
12	Hydrogen Chloride	Isokinetic Sampling ^[5] 1) Absorption Sampling, Ion Chromatographic Method ^[5]
13	Hydrogen Fluoride	2) Isokinetic Sampling, Ion Chromatographic Method ^[5] 1) Absorption Sampling, Ion Chromatographic Method ^[5]
14	Hydrogen Sulfide	2) Isokinetic Sampling, Ion Chromatographic Method ^[5] Absorption Sampling, Iodometric Method ^[5]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
15	Lead	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
16	Manganese	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
17	Mercury	1) Isokinetic Sampling, Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[5] 2) Isokinetic Sampling, Digestion, Cold-Vapor Atomic Fluorescence Spectrometric Method ^[5]
18	Nickel	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
19	Opacity	Ringelmann's Method ^[2]
20	Oxides of Nitrogen	1) Absorption Sampling, Phenoldisulfonic Acid Method ^[5] 2) Absorption Sampling, Alkaline Permanganate/Colorimetric Method ^[5] 3) Instrumental Analyzer Method ^[5]
21	Selenium	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
22	Sulfur Dioxide	1) Absorption Sampling, Barium-Thorin Titrimetric Method ^[5] 2) Instrumental Analyzer Method ^[5]
23	Sulfuric Acid	Isokinetic Sampling, Barium-Thorin Titrimetric Method ^[5]
24	Tellurium	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
25	Tin	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
26	Total Suspended Particulate	1) Isokinetic Sampling, Gravimetric Method ^[5] 2) Paired Train, Isokinetic Sampling, Gravimetric Method ^[5]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
27	Vanadium	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
28	Xylene	Adsorption Sampling, Gas Chromatographic Method ^[5]

สิ่งปฏิกูลหรือวัสดุที่ไม่ใช้แล้ว จำนวน 35 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Aldrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26]
2	Antimony	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]
3	Arsenic	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]
4	Barium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
5	Beryllium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]
6	Cadmium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]
7	Chlordane	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26]
8	Chromium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]
9	Chromium (III)	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method; Waste Extraction, Colorimetric Method; Calculation Method ^[1,6,16,19] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method; Waste Extraction, Colorimetric Method; Calculation Method ^[1,6,17,19] 3) Digestion, Inductively Coupled Plasma Method; Alkaline Digestion, Colorimetric Method; Calculation Method ^[7,8,16,19] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method; Alkaline Digestion, Colorimetric Method; Calculation Method ^[7,8,17,19]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
10	Chromium (VI)	1) Waste Extraction, Colorimetric Method ^[1,6,19] 2) Alkaline Digestion, Colorimetric Method ^[8,19]
11	Cobalt	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]
12	Copper	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]
13	2,4-D	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26]
14	DDD	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26]
15	DDE	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26]
16	DDT	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
17	Dieldrin	2) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 3) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26) 1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^(1,9,26) 2) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 3) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
18	Endrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^(1,9,26) 2) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 3) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
19	Heptachlor	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^(1,9,26) 2) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 3) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
20	Lead	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(1,6,16) 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^(1,6,17) 3) Digestion, Inductively Coupled Plasma Method ^(7,16) 4) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^(7,17)
21	Lindane	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^(1,9,26) 2) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 3) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
22	Mercury	1) Waste Extraction, Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^(1,6,20) 2) Waste Extraction, Digestion, Cold-Vapor Atomic Fluorescence Spectrometric Method ^(1,6,30) 3) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ⁽²⁰⁾ 4) Digestion, Cold-Vapor Atomic Fluorescence Spectrometric Method ⁽³⁰⁾ 5) Thermal Decomposition Amalgamation and Atomic Absorption Spectrometric Method ⁽²¹⁾
23	Methoxychlor	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^(1,9,26) 2) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 3) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
24	Mirex	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^(1,9,26) 2) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 3) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
25	Molybdenum	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(1,6,16) 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^(1,6,17) 3) Digestion, Inductively Coupled Plasma Method ^(7,16) 4) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^(7,17)
26	Nickel	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(1,6,16) 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^(1,6,17) 3) Digestion, Inductively Coupled Plasma Method ^(7,16) 4) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^(7,17)
27	Polychlorinated biphenyls (PCBs) - Aroclor 1016 - Aroclor 1221 - Aroclor 1232 - Aroclor 1242 - Aroclor 1248 - Aroclor 1254 - Aroclor 1260	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^(1,9,26) 2) Soxhlet Extraction, Gas Chromatographic Method ^(10,26) 3) Automated Soxhlet Extraction, Gas Chromatographic Method ^(11,26)

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
28	- 2-Chlorobiphenyl - 2,3-Dichlorobiphenyl - 2,2',5-Trichlorobiphenyl - 2,4',5-Trichlorobiphenyl - 2,2',3,5'-Tetrachlorobiphenyl - 2,2',5,5'-Tetrachlorobiphenyl - 2,3',4,4'-Tetrachlorobiphenyl - 2,2',3,4,5'-Pentachlorobiphenyl - 2,2',4,5,5'-Pentachlorobiphenyl - 2,3,3',4',6-Pentachlorobiphenyl - 2,2',3,4,4',5'-Hexachlorobiphenyl - 2,2',3,4,5,5'-Hexachlorobiphenyl - 2,2',3,5,5',6-Hexachlorobiphenyl - 2,2',4,4',5,5'-Hexachlorobiphenyl - 2,2',3,3',4,4',5-Heptachlorobiphenyl - 2,2',3,4,4',5,5'-Heptachlorobiphenyl - 2,2',3,4,4',5',6-Heptachlorobiphenyl - 2,2',3,4',5,5',6-Heptachlorobiphenyl - 2,2',3,3',4,4',5,5',6-Nonachlorobiphenyl	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26] Electrometric Method ^[23,24]
29	pH	
30	Selenium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
31	Silver	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]
32	Thallium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]
33	Toxaphene	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26]
34	Vanadium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]
35	Zinc	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]

31mg

ดิน จำนวน 125 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Acenaphthene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
2	Acetone	1) Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25) 2) Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method ⁽¹³⁾
3	Aldrin	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
4	Anthracene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
5	Antimony	1) Digestion, Inductively Coupled Plasma Method ^(7,16) 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^(7,17)
6	Arsenic	1) Digestion, Inductively Coupled Plasma Method ^(7,16) 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^(7,17)
7	Atrazine	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
8	Barium	1) Digestion, Inductively Coupled Plasma Method ^(7,16) 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^(7,17)
9	Benz(a)anthracene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
10	Benzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25)

11 Benzo(b)fluoranthene

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
11	Benzo(b)fluoranthene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
12	Benzo(k)fluoranthene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
13	Benzoic acid	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
14	Benzo(a)pyrene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
15	Benzo(g,h,i)perylene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
16	Beryllium	1) Digestion, Inductively Coupled Plasma Method ^(7,16) 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^(7,17)
17	Bis(2-chloroethyl)ether	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
18	Bis(2-ethylhexyl)phthalate	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
19	Bromodichloromethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25)
20	Bromoform	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25)
21	Butanol	Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method ^(13,25)
22	Butyl Benzyl Phthalate	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)

23 Cadmium...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
23	Cadmium	1) Digestion, Inductively Coupled Plasma Method ^[7,16] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]
24	Carbazole	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
25	Carbon Disulfide	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
26	Carbon tetrachloride	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
27	Chlordane	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
28	p-Chloroaniline	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
29	Chlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
30	Chlorodibromomethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
31	Chloroform	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
32	2-Chlorophenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
33	Chromium	1) Digestion, Inductively Coupled Plasma Method ^[7,16] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]
34	Chromium (III)	1) Digestion, Inductively Coupled Plasma Method; Alkaline Digestion, Colorimetric Method; Calculation Method ^[7,8,16,19] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method; Alkaline Digestion, Colorimetric Method; Calculation Method ^[7,8,17,19]
35	Chromium (VI)	Alkaline Digestion, Colorimetric Method ^[8,19]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
36	Chrysene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
37	Cyanide	Extraction, Distillation, Colorimetric Method ^[27,28,29]
38	2,4-D	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
39	DDD	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
40	DDE	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
41	DDT	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
42	Dibenz(a,h)anthracene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
43	Di-n-Butyl Phthalate	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
44	1,2-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
45	1,3-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
46	1,4-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
47	3,3-Dichlorobenzidine	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
48	1,1-Dichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
49	1,2-Dichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25)
50	1,1-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25)
51	cis-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25)
52	trans-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25)
53	2,4-Dichlorophenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
54	1,2-Dichloropropane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25)
55	1,3-Dichloropropane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25)
56	1,3-Dichloropropene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25)
57	Dieldrin	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
58	Diethyl Phthalate	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
59	2,4-Dimethylphenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
60	2,4-Dinitrophenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
61	2,4-Dinitrotoluene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
62	2,6-Dinitrotoluene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
63	Di-n-Octyl Phthalate	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
64	Endosulfan	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
65	Endrin	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
66	Ethylbenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25)
67	Fluoranthene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
68	Fluorene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
69	Heptachlor	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
70	Heptachlor epoxide	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
71	Hexachlorobenzene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
72	Hexachloro-1,3-butadiene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25)
73	n-Hexane	1) Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25) 2) Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method ⁽¹³⁾

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
74	α -HCH	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
75	β -HCH	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
76	γ -HCH	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
77	Hexachlorocyclopentadiene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
78	Hexachloroethane	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
79	Indeno(1,2,3-cd)pyrene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
80	Isophorone	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
81	Lead	1) Digestion, Inductively Coupled Plasma Method ^[7,16] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]
82	Manganese	1) Digestion, Inductively Coupled Plasma Method ^[7,16] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]
83	Mercury	1) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[20] 2) Thermal Decomposition, Amalgamation, and Atomic Absorption Spectrophotometry ^[21] 3) Digestion, Cold-Vapor Atomic Fluorescence Spectrometric Method ^[30]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
84	Methanol	1) Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25] 2) Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
85	Methoxychlor	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
86	Methyl Bromide	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
87	Methylene Chloride	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
88	2-methylphenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
89	2-Methylnaphthalene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
90	Methyl tert-Butyl Ether	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
91	Naphthalene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
92	Nickel	1) Digestion, Inductively Coupled Plasma Method ^[7,16] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]
93	Nitrobenzene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
94	N-Nitrosodiphenylamine	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
95	N-Nitrosodi-n-propylamine	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
96	Polychlorinated biphenyls (PCBs) - Aroclor 1016 - Aroclor 1221 - Aroclor 1232 - Aroclor 1242 - Aroclor 1248 - Aroclor 1254 - Aroclor 1260 - 2-Chlorobiphenyl - 2,2',3,5'-Tetrachlorobiphenyl - 2,2',5,5'-Tetrachlorobiphenyl - 2,3',4,4'-Tetrachlorobiphenyl - 2,2',3,4,5'-Pentachlorobiphenyl - 2,2',4,5,5'-Pentachlorobiphenyl - 2,3,3',4',6-Pentachlorobiphenyl - 2,2',3,4,4',5'-Hexachlorobiphenyl - 2,2',3,4,5,5'-Hexachlorobiphenyl - 2,2',3,5,5',6-Hexachlorobiphenyl - 2,2',4,4',5,5'-Hexachlorobiphenyl - 2,2',3,3',4,4',5-Heptachlorobiphenyl - 2,2',3,4,4',5,5'-Heptachlorobiphenyl - 2,2',3,4,4',5,6-Heptachlorobiphenyl - 2,2',3,4',5,5',6-Heptachlorobiphenyl - 2,2',3,3',4,4',5,5',6-Nonachlorobiphenyl	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
97	Pentachlorophenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
98	Phenanthrene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
99	Phenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
100	Pyrene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
101	Selenium	1) Digestion, Inductively Coupled Plasma Method ^[7,16] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]
102	Silver	1) Digestion, Inductively Coupled Plasma Method ^[7,16] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]
103	Styrene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
104	1,1,2,2-Tetrachloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
105	Tetrachloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
106	Toluene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
107	Toxaphene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
108	TPH (C ₅ -C ₈)	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
109	TPH (C ₈ -C ₁₆)	1) Automate Extraction, Gas Chromatographic Method ^[11,22] 2) Solvent Extraction, Gas Chromatographic Method ^[11,22] 3) Ultrasonic Extraction, Gas Chromatographic Method ^[22,31]
110	TPH (C ₁₆ - C ₃₅)	1) Automate Extraction, Gas Chromatographic Method ^[11,22] 2) Solvent Extraction, Gas Chromatographic Method ^[11,22] 3) Ultrasonic Extraction, Gas Chromatographic Method ^[22,31]
111	1,2,4-Trichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
112	1,1,1-Trichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
113	1,1,2-Trichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
114	Trichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
115	2,4,5-Trichlorophenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
116	2,4,6-Trichlorophenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(10,26) 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^(11,26)
117	1,3,5-Trimethylbenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25)
118	Vanadium	1) Digestion, Inductively Coupled Plasma Method ^(7,16) 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^(7,17)
119	Vinyl Acetate	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25)
120	Vinyl Chloride	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25)
121	m-Xylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25)
122	o-Xylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25)
123	p-Xylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25)
124	Xylene (Total)	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^(15,25)
125	Zinc	1) Digestion, Inductively Coupled Plasma Method ^(7,16) 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^(7,17)

เอกสารอ้างอิง

- กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2566. เรื่อง การจัดการสิ่งปฏิกูลหรือวัสดุที่ไม่ใช้แล้ว. ราชกิจจานุเบกษา. 31 พฤษภาคม 2566. เล่มที่ 140 ตอนพิเศษ 126 ง.
- กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2549. เรื่อง กำหนดค่าปริมาณเข้มข้นที่เจือปนในอากาศที่ระบายออกจากปล่องของหม้อน้ำโรงสีข้าวที่ใช้แก๊สเป็นเชื้อเพลิง. ราชกิจจานุเบกษา. 4 ธันวาคม 2549. เล่มที่ 123 ตอนพิเศษ 125ง.
- สมาคมวิศวกรรมสิ่งแวดล้อมแห่งประเทศไทย. คู่มือวิเคราะห์น้ำเสีย. พิมพ์ครั้งที่ 4. กรุงเทพฯ: เรือนแก้วการพิมพ์, 2547.
- APHA, AWWA, WEF. Standard Methods for the Examination of Water and Wastewater. 24th ed. Washington, DC: APHA, 2023.

5. United States...

- United States Environmental Protection Agency. Standards of Performance for New Stationary Sources. 40 CFR 60. Appendix A, 2023.
- United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. SW-846, 2014.
- United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Acid Digestion of Sludges and Sediments and Soils. SW-846 Method 3050B, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Alkaline Digestion for Hexavalent Chromium. SW-846 Method 3060A, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Separatory Funnel Liquid-Liquid Extraction. SW-846 Method 3510C, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Soxhlet Extraction. SW-846 Method 3540C, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Automated Soxhlet Extraction. SW-846 Method 3541, 1994.
- United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Microscale Solvent Extraction (MSE). SW-846 Method 3570, 2002.
- United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Volatile Organic Compounds (VOCs) in Various Sample Matrices Using Equilibrium Headspace Analysis. SW-846 Method 5021A, 2014.
- United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Purge-and-Trap for Aqueous Samples. SW-846 Method 5030B, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Closed-System Purge-and-Trap and Extraction for Volatile Organics in Soil and Waste Samples. SW-846 Method 5035, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Inductively Coupled Plasma-Atomic Emission Spectrometry. SW-846 Method 6010B, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Inductively Coupled Plasma-Mass Spectrometry. SW-846 Method 6020A, 2007.
- United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Antimony and Arsenic (Atomic Absorption, Borohydride Reduction). SW-846 Method 7062, 1994. เพิ่มใหม่
- United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Chromium, Hexavalent (Colorimetric). SW-846 Method 7196A, 1992.

20. United States...

20. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Mercury In Solid or Semisolid Waste (Manual Cold-Vapor Technique). SW-846 Method 7471B, 2007.

21. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Mercury in Solids and Solutions by Thermal Decomposition, Amalgamation, and Atomic Absorption Spectrophotometry. SW-846 Method 7473, 2007.

22. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Nonhalogenated Organics by Gas Chromatography. SW-846 Method 8015C, 2007.

23. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. pH Electrometric Measurement. SW-846 Method 9040C, 2004.

24. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Soil and Waste pH. SW-846 Method 9045D, 2004.

25. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Volatile Organic Compounds by Gas Chromatography/Mass Spectrometry (GC/MS). SW-846 Method 8260D, 2018.

26. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Semivolatile Organic Compounds by Gas Chromatography/Mass Spectrometry (GC/MS). SW-846 Method 8270E, 2018.

27. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Total and Amenable Cyanide: Distillation SW-846 Method 9010B, 1996.

28. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Cyanide Extraction Procedure for Solids and Oil. SW-846 Method 9013A, 1996.

29. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Cyanide in Waters and Extracts Using Titrimetric and Manual Spectrophotometric Procedures. SW-846 Method 9014, 2014.

30. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Mercury in Sediment and Tissue Samples by Atomic Fluorescence Spectrometry. SW-846 Method 7474, 2007.

31. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Ultrasonic Extraction. SW-846 Method 3550C, 2007.

Ampl



ที่ อก ๐๓๑๐(๑)/ ๔ ๑ ๒ ๑

กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๕๐๐

๒๕ เมษายน ๒๕๖๗

เรื่อง เปลี่ยนแปลงบุคลากรของห้องปฏิบัติการวิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๒๙ มีนาคม ๒๕๖๗

ตามคำขอที่อ้างถึง บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด ห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ว-๒๐๔ สถานที่ตั้งเลขที่ ๑๐๔ ซอยพัฒนาการ ๔๐ ถนนพัฒนาการ แขวงพัฒนาการ เขตสวนหลวง กรุงเทพมหานคร ขอเปลี่ยนแปลงบุคลากร ความละเอียดแจ้งแล้ว นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว มีความเห็นดังนี้

๑. ให้ยกเลิกเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๓ ราย

๑) นางสาวพรรณธิดา พุ่มคง	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๕
๒) นายกำชัย สุทธระ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๒๑
๓) นางสาวศุภรดา ปันมยุรา	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๓๘

๒. ให้เพิ่มเจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกชน จำนวน ๑๒ ราย

๑) นางสาวฐานิดา กลิ่นเขียว	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๒
๒) นางสาวกัญญ์กิสสร สายคำ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๓
๓) นางสาวณัฐนันท์ กันทะวงศ์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๔
๔) นายอำนาจ วงษาเคน	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๕
๕) นายกฤษณพล ปัญญาวงศ์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๖
๖) นายณชากร หารรักษา	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๗
๗) นายวัชรินทร์ ผ่องสามสวน	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๘
๘) นายณัฐพงศ์ โสภา	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๙
๙) นายศักรินทร์ ปานเพ็ง	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๙๐
๑๐) นายณัฐพล ชุ่มชื่น	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๙๑
๑๑) นายธนา สุภาพันธ์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๙๒
๑๒) นายณรรธ แก้วพงษ์ชา	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๙๓

อนึ่ง หนังสือฉบับนี้...

อนึ่ง หนังสือฉบับนี้จะหมดอายุพร้อมหนังสือต่ออายุรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
ในวันที่ ๒ กันยายน ๒๕๖๙

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

พ

(นายพรยศ กลั่นกรอง)

รองอธิบดี ปฏิบัติราชการแทน
อธิบดีกรมโรงงานอุตสาหกรรม

กองวิจัยและเตือนภัยมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕

โทรสาร ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๕๙

ไปรษณีย์อิเล็กทรอนิกส์ saraban@div.mail.go.th

ที่ อก ๐๓๒๒/๑๓๖๕๙



กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๒๕ ก.ย. ๒๕๖๖

เรื่อง ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

เรียน กรรมการผู้จัดการ บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๑๔ กรกฎาคม ๒๕๖๖

สิ่งที่ส่งมาด้วย เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด จำนวน ๓ แผ่น

ตามหนังสือที่อ้างถึง บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด ขอต่ออายุ
หนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ว-๒๖๗ สถานที่ตั้ง เลขที่ ๑๑๔/๑ หมู่ที่ ๘
ถนนกาญจนาภิเษก ตำบลบ้านพรุ อำเภอหาดใหญ่ จังหวัดสงขลา ต่อกรมโรงงานอุตสาหกรรม นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย)
จำกัด ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน โดยมีองค์ประกอบดังนี้

ก. ผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์

นางสาวกนิษฐา เหมประสาทร

ทะเบียนเลขที่ ว-๒๖๗-ค-๐๐๐๑

ข. เจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์

๑) นางสาวอินทิรา คงประยูร

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๐๑

๒) นางสาวอมรรัตน์ เพชรประดับ

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๐๒

๓) นายทักษิณ อินโดรม

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๐๓

๔) นางสาวอนันดา บุญเพชร

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๐๔

๕) นางสาวสุทธิรักษ์ ทิพย์รัตน์

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๐๕

๖) นางสาวนริสา นฤมิตร

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๐๖

๗) นายวุฒิชัย ทวยเจริญ

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๐๗

๘) นายยงศิลป์ รั้งชัย

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๐๘

๙) นายอภิวัฒน์ ฉันทะ

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๐๙

๑๐) นายศิริชัย เกสียงเกิด

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๑๐

๑๑) นายสมศักดิ์ จันทรงค์

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๑๑

๑๒) นางสาวพิชญา ศุภรานนท์

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๑๒

๑๓) นายปัญญา เกียรติพิริรักษ์

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๑๓

๑๔) นางสาวศศิณิภา รอดทองอ่อน

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๑๔

๑๕) นางสาวชุติมา สุขสวัสดิ์

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๑๕

๑๖) นางสาวจันทิมา คงทน

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๑๖

๑๗) นางสาวกุลวดี เรืองประพันธ์

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๑๗

๑๘) นางสาวอาทิตย์ยา เมืองแก้ว

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๑๘

๑๙) นางสาวกวิณา ฉุนย่อง

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๑๙

ค. ขอบข่ายสารมลพิษที่ได้รับขึ้นทะเบียนให้วิเคราะห์ในน้ำเสียและอากาศเสีย ตามสิ่งที่ส่งมาด้วย

หนังสือฉบับนี้...



หนังสือฉบับนี้จะหมดอายุในวันที่ ๑๖ สิงหาคม ๒๕๖๙ หากประสงค์จะต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ให้อื่นคำขอต่ออายุพร้อมเอกสารประกอบคำขอต่อกรมโรงงานอุตสาหกรรมภายใน ๓๐ วัน ก่อนวันสิ้นอายุของหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ทั้งนี้ สามารถยื่นคำขอผ่านระบบอิเล็กทรอนิกส์ได้ที่หน้าเว็บไซต์กรมโรงงานอุตสาหกรรม

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

น

(นายณเรศวร์ ตริยงค์)

ผู้อำนวยการศูนย์วิจัยและเตือนภัยมลพิษโรงงานภาคใต้
ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม

ศูนย์วิจัยและเตือนภัยมลพิษโรงงานภาคใต้
โทร. ๐ ๗๔๓๒ ๕๐๒๙, ๐ ๗๔๘๙ ๐๖๓๔ ต่อ ๕๒๐๑
ไปรษณีย์อิเล็กทรอนิกส์ sirw@diw.mail.go.th

เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
บริษัท เอแอลเอส แลบลอราทอรี กรุ๊ป (ประเทศไทย) จำกัด เลขทะเบียน ว-๒๖๗
ที่ อก ๐๓๒๒/๑๗๖๕๔/ ลงวันที่ ๒๕ ก.ย. ๒๕๖๖

ขอขยายสารมลพิษที่ได้รับการขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๓๗ รายการ
น้ำเสีย จำนวน 25 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Arsenic	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1]
2	Barium	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1]
3	Biochemical Oxygen Demand	5-Day BOD Test, Azide Modification Method ^[1] 5-Day BOD Test, Membrane Electrode Method ^[1]
4	Cadmium	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1]
5	Chemical Oxygen Demand	Closed Reflux, Colorimetric Method ^[1] Closed Reflux, Titrimetric Method ^[1]
6	Chromium	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1]
7	Color	ADMI Weighted-Ordinate Spectrophotometric Method ^[1]
8	Copper	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1]
9	Formaldehyde	Distillation, Colorimetric Method ^[2]
10	Free Chlorine	DPD Ferrous Titrimetric Method ^[1]
11	Hexavalent Chromium	Filtration, Colorimetric Method ^[1]
12	Lead	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1]
13	Manganese	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1]
14	Mercury	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1]
15	Nickel	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1]
16	Oil & Grease	Liquid-Liquid, Partition-Gravimetric Method ^[1]

นางสาว รติณสุภา
(นางสาวบุษยา รติณสุภา)
นักวิทยาศาสตร์ชำนาญการ

17 pH...

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
17	pH	Electrometric Method ^[1]
18	Phenol	Distillation, Direct Photometric Method ^[1]
19	Selenium	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[1]
20	Sulfide	ZnS Precipitation, Iodometric Method ^[1]
21	Temperature	Laboratory and Field Methods ^[1]
22	Total Dissolved Solids	Dried at 180 °C ^[1]
23	Total Suspended Solids	Dried at 103-105 °C ^[1]
24	Trivalent Chromium	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[1]
25	Zinc	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[1]

อากาศเสีย จำนวน 12 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Antimony	Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[3]
2	Arsenic	Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[3]
3	Carbon Monoxide	Sampling Bag Non-Dispersive Infrared Method ^[3]
4	Copper	Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[3]
5	Dioxins	Isokinetic Sampling, Analysis by ISO/IEC 17025 Accredited Laboratory or Analysis by Department of Industrial Works Registered Laboratory ^[3]
6	Hydrogen Sulfide	Absorption, Iodometric Method ^[3]
7	Lead	Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[3]
8	Opacity	Ringelmann's Method ^[4]
9	Oxides of Nitrogen	Absorption Sampling, Phenoldisulfonic acid Method ^[3]
10	Sulfur Dioxide	Absorption Sampling, Barium-Thorin Titrimetric Method ^[3]
11	Sulfuric acid	Isokinetic Sampling, Barium-Thorin Titrimetric Method ^[3]
12	Total Suspended Particulate	Isokinetic Sampling, Gravimetric Method ^[3]

บุษยา รัตนสุภา
(นางสาวบุษยา รัตนสุภา)
นักวิทยาศาสตร์ชำนาญการ

เอกสารอ้างอิง....

เอกสารอ้างอิง

1. APHA, AWWA, WEF. Standard Methods for the Examination of Water and Wastewater. 24th ed. Washington, DC: APHA, 2023.
2. สมาคมวิศวกรรมสิ่งแวดล้อมแห่งประเทศไทย. คู่มือวิเคราะห์น้ำเสีย. พิมพ์ครั้งที่ 4. กรุงเทพฯ: เรือนแก้วการพิมพ์, 2547.
3. United States Environmental Protection Agency. Standards of Performance for New Stationary Sources. 40 CFR 60. Appendix A, 2020.
4. กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2549. เรื่อง กำหนดค่าปริมาณเขม่าควันที่เจือปนในอากาศที่ระบายออกจากปล่องของหม้อน้ำโรงสีข้าวที่ใช้แก๊สเป็นเชื้อเพลิง. ราชกิจจานุเบกษา. 4 ธันวาคม 2549. เล่มที่ 123 ตอนพิเศษ 125ง.

บุษยา รัตนสุภา
(นางสาวบุษยา รัตนสุภา)
นักวิทยาศาสตร์ชำนาญการ



right solutions.
right partner.

บริษัท เอแอลเอส แล็บอราทอรี กรุ๊ป (ประเทศไทย) จำกัด (สำนักงานใหญ่)
104 ซอยพัฒนาการ 40 ถนนพัฒนาการ
แขวงพัฒนาการ เขตสวนหลวง กรุงเทพฯ 10250



ติดต่อเรา

