

ภาคผนวก ช
เอกสารสอบเทียบเครื่องมือ



List of Instruments Certification for Water Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration	Remark
Water									
1	pH Meter	pH	Horiba	LAQUA-PH210 HA0C0025	Technology Promotion Association (Thailand-Japan)	24CH319	14 Mar 24	13 Mar 25	-



Cert.No.: 24CH319
Page: 1 of 3

Certificate of Calibration

Equipment : pH Meter
Manufacturer : HANNA
Model : I-AQUA-PH210
Serial No. : HA000025
ID No. : UA-EFM-117256Q/EFM-pH137853
Condition As-Received: Used Item
Received Date : 12 March 2024
Calibration Date : 14 March 2024
Reference : 2403 03B3W5C 1
Submitted by : United Analyst and Engineering Consultant Co., Ltd.
381 Udomsuk 41, Sukhumvit Road,
Bangnaek, Phraekhanong, Bangkok 10760
Ambient Temperature : (25 ± 2.5) °C
Relative Humidity : (50 ± 15) %
Calibration Procedure :
- CP-C1-5 by direct measurement with DC voltage
standard and direct measurement with
certified reference materials (CRM)
- CP-C1-8 by comparison with temperature standard

Calibrated by : Warakorn Lungsakulakul

Approved by :

() Pankhappa Tamsayaku
() Unnopphol Harnachai
(✓) Sathip Boonngima

Issue Date : 15 March 2024

The Uncertainties are for a confidence probability of approximately 95 %

This certificate is the property of the user for whom it is issued and is not to be
reproduced or used for any other purpose without the written consent of the
issuing organization.

เอกสารไม่ควบคุม

A 0064529



Cert.No.: 24CH319
Page: 2 of 3

Condition of this calibration result

Reference Standard Instrument

Instrument	Serial No.	ID No.	Cert. No.	Exp. Date
1) Document Process Calibrator	54030048	150407118	23CE002	27 Aug 2024
2) Ref. Standard Thermometer	4982054	11040044	231908	26 July 2024

This calibration is traceable to the International System of Unit maintained through
Technology Promotion Association (Thailand-Japan)

2. Certified Reference Materials : The measurement results are traceable to 5 through CPA chem Ltd.,
AMS-ASQ National Accreditation Board, Accredited No. AS-1928

Buffer Solution	Manufacturer	Lot No.	Exp. date
pH 4.008	CPA chem	940102	27 Nov 2025
pH 6.865	CPA chem	940104	02 Nov 2024
pH 9.997	CPA chem	940106	02 Nov 2024

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration Results

Function : mV Measurement

Performing standard curve by Document Process Calibrator at pH (4.7)(7,10)

Unit Under Calibration	Nominal Value	Standard Voltage Input	Actual Reading		Uncertainty of Measurement (± mV)	Coverage factor k
	pH	mV	mV	pH		
pH Meter S/N.: HA000025	4.00	177.48	177.5	4.01	0.028	2.00
	7.00	0.00	0.0	7.03	0.058	2.00
	7.00	0.00	0.0	7.03	0.058	2.00
	10.00	-177.48	-177.5	10.04	0.058	2.00

เอกสารไม่ควบคุม

a 1206341



Cert.No.: 24CH219
Page: 3 of 3

Calibration Results

Function : pH Measurement

Performing three buffers standard curve by using buffer nominal pH (4.7)(7,10)

Unit Under Calibration	Standard pH Buffer Solution	Actual pH Reading	Actual mV Reading (mV)	Uncertainty of pH Measurement (±)	Coverage factor k
pH Electrode S/N.: -	4.008	4.01	145.4	0.0021	2.27
	6.866	7.00	-25.1	0.0035	2.00
	6.866	7.02	24.3	0.014	2.00
	9.997	9.01	-149.5	0.0035	2.00

Function : Temperature Measurement

(*) Without adjustment

This equipment was connected with Temperature Probe

- Model :
- Serial No. :
Dimension of probe
- Length : 120 mm.
- Diameter : 16 mm.
- Insertion Depth : 60 mm.

Calibration Point (°C)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty of measurement (± °C)	Coverage factor k
25.0	25.001	25.0	-0.021	0.13	2.00
20.0	20.001	20.0	-0.021	0.13	2.00
25.0	25.002	25.0	-0.022	0.12	2.00

Remark : *UUC = Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-o-o-

เอกสารไม่ควบคุม

a 1206342

List of Instruments Certification for Environmental Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*	Remark
Equipment for Soil Quality Analysis									
1	pH Meter	ความเป็นกรดด่าง (pH)	Mettler-Toledo	Seven Easy S20 / 1231155210	National Food Institute, Ministry of Industry, Thailand	2401718-001-01	11 Mar 24	10 Mar 25	-
2	pH Meter		Mettler-Toledo	Seven Easy S20 / 1230525212	DKSH (Thailand) Ltd.	C07240167	9 Apr 24	8 Apr 25	-
3	Analytical Balance (Readability 0.1 mg)	สารกลุ่มปิโตรเลียมไฮโดรคาร์บอนทั้งหมด (TPH)	Mettler-Toledo	AB-204S/FACT / 1129361010	Technology Promotion Association (Thailand-Japan)	24MM292	11 May 24	10 May 25	-
4	Atomic Absorption Spectrophotometer (AAS)	กลุ่มโลหะหนัก (As, Cd, Pb, Hg, Ba)	Agilent Technologies	System ID:G8432A AA240FS / MY13160001	Thailand Institute of Scientific and Technological Research(TISTR)	MTC.ACL.No 358/67	11 Mar 24	10 Mar 25	-
5	Inductively Coupled Plasma (ICP)		Agilent Technologies	System ID:G8015A G8015AA / MY18030001	Agilent Technologies (Thailand) Co.,Ltd.	Preventive Maintenance Checklist	4 Nov 24	3 Nov 25	-

Due Date of Calibration* : Based on the annual calibration plan. At least 1 time per year.

Calibration Certificate

Certificate No.: 2401718-001-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road,
Bangchack, Prakanong, Bangkok 10260

Page 1 of 5

Equipment: pH Meter
Manufacturer: METTLER TOLEDO
Model: SevenEasy pH
Serial No.: 1231155210
ID No.: UAE.WAT.010/2553
Order No.: 2401718
Operation No.: 2401718-001
Date of Receipt: 27 February 2024
Date of Calibration: 11 March 2024

Calibrated by Mr.Manas Somsak Specialist
Approved by (Mr.Pheraphat Tuanjit)
Manager, Division of Calibration Laboratory
Responsible for the Technical Management Team
Date of Issue: 12 March 2024

The uncertainties are for a confidence probability of approximately 95%.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65



Calibration Report

Certificate No.: 2401718-001-01
Equipment: pH Meter
Resolution: 0.01 pH : 1 mV
Manufacturer: METTLER TOLEDO
Model: SevenEasy pH
Serial No.: 1231155210
Type: Bench top
ID No.: UAE.WAT.010/2553

Date of Calibration: 11 March 2024 Page 3 of 5

1. Calibration of pH Meter (Manual Temperature Compensation at 25 °C)

Nominal pH	DC Voltage Standard (mV)	Average Indicator Reading		Uncertainty (±mV)	Coverage Factor (k)
		mV	pH		
0	414.121	414	0.00	0.58	2.00
2	296.814	296	2.00	0.58	2.00
4	177.464	178	4.00	0.58	2.00
6	59.160	59	6.00	0.58	2.00
7	0.001	0	7.00	0.58	2.00
8	-59.159	-59	8.00	0.58	2.00
10	-177.461	-177	10.00	0.58	2.00
12	-296.811	-296	12.00	0.58	2.00
14	-414.118	-414	14.00	0.58	2.00

2. Calibration of pH Meter with Electrode (Manual Temperature Compensation at 25 °C)

Equipment: pH Electrode Type: Combined Electrode
Manufacturer: METTLER TOLEDO Model: InLab Solids
Serial No.: 3065701 ID No.: N/A

Performance of Electrode system (Three-Point Calibration at pH 4, 7 and 10)

Certified Value (25 °C (pH))	Average Indicator Reading		Relative Slope (%)	Uncertainty (± pH)	Coverage Factor (k)
	pH	mV			
4.000	4.01	188	-	0.0071	2.00
7.001	7.00	13	98.9	0.0086	2.00
10.010	10.01	-160	97.2	0.0085	2.00
6.865	6.87	21	-	0.0074	2.00

F-CS-012 Revision: 01 Date: 20-04-65



Calibration Report

Certificate No.: 2401718-001-01
Equipment: pH Meter
Resolution: 0.01 pH : 1 mV
Manufacturer: METTLER TOLEDO
Model: SevenEasy pH
Serial No.: 1231155210
Type: Bench top
ID No.: UAE.WAT.010/2553

Date of Calibration: 11 March 2024 Page 2 of 5

Location: Chemical Calibration Laboratory, National Food Institute
Environment Condition: Ambient Temperature: (23.4 ± 1.5) °C Relative Humidity: (51 ± 3) %
Condition of Equipment: Good Condition
Condition of this Results of Calibration

1. Calibration Method: W-CC-002 : In house method based on direct measurement by using standard voltage calibrator and certified reference material (CRM)

2. Reference Standards / Certified Reference Material

Instruments	Serial / ID No.	Manufacturer	Certificate No.	Due Date	
2.1 DC Voltage Calibrator	2709007	Fuke	Z0E2003	14 June 2024	
2.2 Digital Thermometer	2709007	Fuke	CC 660570-01	30 October 2024	
2.3 Thermo-Hygro Meter	NFLBTH 014023	teiso	CC 660353-01	3 April 2024	
Certified Reference Material		Lot No.	Manufacturer	Ref N	Expiry Date
2.4 pH buffer 4.008 (Primary pH buffer Solution)	868842	CPAchem	PH216 L5		13 April 2025
2.5 pH buffer 6.865 (Primary pH buffer Solution)	868843	CPAchem	PH217 L5		13 April 2025
2.6 pH buffer 10.01 (Primary pH buffer Solution)	868844	CPAchem	PH220 L5		13 April 2024
2.7 pH buffer 7.00 (Standard pH buffer Solution)	033109	HACH LANGE GmbH	S11M004		16 October 2025

3. This certification is traceable to The International System of Unit (SI Unit)

- 3.1 Instruments Ng.2.1 through NSC-TISI-TS 17025 Laboratory Accreditation of Calibration No.0008
- 3.2 Instruments Ng.2.2 and 2.3 through NSC-TISI-TS 17025 Laboratory Accreditation of Calibration No.0061
- 3.3 Certified Reference Material Ng.2.4 to 2.6 traceable to Primary measurement method- Hamed cell using calibrated thermometer, barometer, and nanovisimeter. The Standard Solution preparation and certified by CPAchem Ltd is accredited to ISO 17034 and ISO/IEC 17025
- 3.4 Certified Reference Material Ng.2.7 traceable to PTB Certificate Nr. PTB-PH04-563/2009423 and Certificate Nr. PTB-PH08-555/3062022 (PTB: Physikalisch-Technische Bundesanstalt, Braunschweig, Germany)

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

F-CS-012 Revision: 01 Date: 20-04-65



Calibration Report

Certificate No.: 2401718-001-01
Equipment: Digital Thermometer with RTD (pH Meter)
Resolution: 0.1 °C Model: SevenEasy pH
Serial No.: 1231155210 ID No.: UAE.WAT.010/2553
Manufacturer: METTLER TOLEDO

Date of Calibration: 11 March 2024 Page 4 of 5

Location: Chemical Calibration Laboratory, National Food Institute

Environment Condition: Ambient Temperature 23 °C ± 1 °C
Relative Humidity 51 % ± 2 %

Condition of this results of Calibration:

- 1. Calibration Method : - In house method: W-TE-025 by comparison with standard thermometer
- The Calibration is determined by comparing with a known temperature from a standard resistance thermometer.
- The temperature scale in use at this laboratory is the International Temperature scale of 1990 (ITS-90).

2. Reference Standard Instrument :

Instrument	Model	Serial No.	Certificate No.	Due Date	Through
HANDHELD THERMOMETER	1523	2118154	PSL-T 087766	06-Jun-24	TSTR
Platinum Resistance Thermometer (PRT)	5627A	877332			

Support Equipment : -Low Temperature Bath (ISOCAL-6), Model: Europa-6 Plus Basic, SN: 341582/2

3. This certificate is traceable to International System of Units (SI Units).

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

6. Condition of Calibrated Item :

7. Result of Calibration : Good
☒ Without adjustment ☐ After adjustment

F-CS-012 Revision: 01 Date: 20-04-65



Calibration Report

Certificate No.: 2401718-001-01
Equipment: Digital Thermometer with RTD (pH Meter)
Resolution: 0.1 °C Model: SevenEasy pH
Serial No.: 1231156210 ID No.: UAE.WAT.010/0553
Manufacturer: METTLER TOLEDO
Date of Calibration: 11 March 2024 Page 5 of 5

Calibration point: 15.0, 25.0 and 35.0 °C

Calibration result:
- The probe was immersed in liquid bath or dry bath to a minimum depth of 100 mm.
- Description of probe, model: N/A S/N: N/A
Dimension of probe: Diameter 4 mm, Length 120 mm.
Sheath material: Stainless Steel

UUC Reading (°C)	Standard Temperature (°C)	Correction Value (°C)	Uncertainty ± (°C)
15.1	14.998	0.1	0.099
25.1	24.998	0.1	0.099
35.1	34.997	0.1	0.099

Note:
- UUC: Unit Under Calibration

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor $k=2$, providing a level of confidence of approximately 95 %.

End

F-CS-012 Revision: 01 Date: 20-04-85



เอกสารไม่ควบคุม

nfi.co.th



Certificate No.: C07240167 Page 2 of 3

Calibration Results:

pH Scale

Input (mV)	pH Meter Reading			Uncertainty of Measurement (mV)	Coverage Factor (k)
	(mV)	Error (mV)	(pH)		
414.12	414	-0.12	0.00	0.58	2.00
354.96	355	0.04	1.00	0.58	2.00
295.8	296	0.20	2.00	0.58	2.00
236.64	237	0.36	3.00	0.58	2.00
177.48	178	0.52	4.00	0.58	2.00
118.32	118	-0.32	5.00	0.58	2.00
59.16	59	-0.16	6.00	0.58	2.00
0	0	0.00	7.00	0.58	2.00
-59.16	-59	0.16	8.00	0.58	2.00
-118.32	-118	0.32	9.00	0.58	2.00
-177.48	-177	0.48	10.00	0.58	2.00
-236.64	-236	0.64	11.00	0.58	2.00
-295.8	-296	-0.20	12.00	0.58	2.00
-354.96	-355	-0.04	13.00	0.58	2.00
-414.12	-414	0.12	14.00	0.58	2.00

Certificate of Calibration

Equipment: pH METER Certificate No.: C07240167
Model: SevenEasy Issued Date: 9 April 2024
Serial No. (or ID.): 1230525212 (UAE.WAS.003/2553) Job No.: WO-00024208
Manufacturer: METTLER TOLEDO Page: 1 of 3
Electrode Serial No.: 1156883 Model: InLab Solids Brand: METTLER TOLEDO
Condition: In Condition

Customer: United Analyst and Engineering Consultant Company Limited
3 Soi Udomsuk 41 Sukhumvit Road,
Bangchak, Prakanong, Bangkok 10260 Thailand

Environment Condition: Temperature 23 °C ± 2 °C
Humidity 50 %RH ± 15 %RH

Calibration Place: Environment Laboratory, DKSH Technology Limited,
2533 Sukhumvit Road, Bangchak,
Phrakanong, Bangkok 10260 Thailand

Calibration By: Miss.Orawan Khlaiphloi

Calibration Date: 9 April 2024

The Method used: In house method, CAL-WI-58, base on ASTM E 70-07

Traceability: This certificate is traceable to SI Units, Sample Test is assured through primary measurement method Harned cell, through CPAchem Ltd. (ISO/IEC 17034) Certificate No. 938377, 931985, 931984 And pH Scale traceable to the SI Units maintained by National Institute of Metrology (NIMT), Thailand through Industrial Foundation Electrical and Electronics Institute Certificate No. CA20230350EA

Orawan
(Miss Orawan Khlaiphloi)
Person in charge

Nitinun
(Mr. Nitinun Srihawan)
Authorized signatory

This certificate is issued the units of measurement according to the International System of Units (SI). It provides traceability of measurement to international or national standard or other recognized national standard laboratories.
The measurement uncertainty stated is the expanded uncertainty which is obtained from the standard uncertainty multiplied by the coverage factor ($k=2$) to provide a level of confidence of approximately 95%. It is determined in accordance with the Guide to Expression of Uncertainty in Measurement (GUM).
These results may be affected by deviations from specified conditions. The results relate only to the items tested, calibrated or sampled. The report shall not be reproduced except in full without approval of DKSH Technology Limited.

DKSH Technology Limited
2533 Sukhumvit Road, Bangkok, Prakanong, Bangkok 10260
Phone: +66 2639 7000 Email: info.calibration@dksh.com Website: www.dksh.com/scientific-thailand

Delivering Growth - in Asia and Beyond.

เอกสารไม่ควบคุม

CAL-FM-C07-14; 9 Apr 2024



Certificate No.: C07240167 Page 3 of 3

Practical slope and zero point*

The three-point calibration using three standard buffer solutions; pH 4.008, pH 6.985 and pH 9.997

-During calibration, display of pH meter reading; pH 4.00, pH 7.00 and pH 10.01

The practical slope of the pH electrode; 57.01 (mV/pH), 96.37%

The zero point of the pH electrode; 6.88 (pH)

Sample Test Results

Standard Buffer Solution (pH)	Unit Under Calibration (pH)	Difference (pH)	Uncertainty of Measurement (pH)	Coverage Factor (k)
4.008	3.99	-0.018	0.0070	2.00
6.985	7.00	0.015	0.0091	2.00
9.997	10.02	0.023	0.0074	2.00

* Calibration Marked * Not TISI Accredited * in this Certificate have been included for completeness.

The End of Certificate



Certificate of Calibration

Certificate No.: C15240373

Page: 2 of 2

Equipment: Digital Thermometer with Probe Certificate No.: C15240373
Model: SevenEasy pH Issued Date: 09 April 2024
Serial No.: 1230525212 Job No.: WO-00024208
Manufacturer: METTLER TOLEDO Page: 1 of 2
ID No.: UAE.WAS.003/2553 Condition: In Condition

Customer: United Analyst and Engineering Consultant Company Limited
3 Soi Udomsuk 41 Sukhumvit Road,
Bangkok, Prakanong, Bangkok 10260 Thailand

Environment Condition: Temperature: 22 °C ± 3 °C
Humidity: 50 %RH ± 20 %RH
Voltage: 220 VAC ± 10 %

Calibration Place: Thermo-Hygro Laboratory, DKSH Technology Limited,
2533 Sukhumvit Road, Bangkok,
Phrakhanong, Bangkok 10260 Thailand

Calibration By: Mr. Nateekarn Mitjit
Calibration Date: 09 April 2024
The Method used: In house method, CAL-WI-19, by comparison with standard thermometer
Traceability: This certificate is traceable to the International System of Unit maintained by
Quality Reborn Co.,Ltd. (QR) Certificate No. QR23-1073

(Mr. Nateekarn Mitjit)

Person in charge

(Mr. Pramote Ramrong)

Authorized signatory

This certificate is issued the units of measurement according to the International System of Units (SI). It provides traceability of measurement to international or national standard or other recognized national standard laboratories.
The measurement uncertainty stated is the expanded uncertainty which is obtained from the standard uncertainty multiplied by the coverage factor (k=2) to provide a level of confidence of approximately 95%. It is determined in accordance with the Guide to Expression of Uncertainty in Measurement (GUM).
These results may be affected by deviations from specified conditions. The results relate only to the items tested, calibrated or sampled. The report shall not be reproduced except in full without approval of DKSH Technology Limited.

บริษัท ดีเคเอส อีเซีย จำกัด
DKSH Technology Limited
2533 ถนนสุขุมวิท แขวงคลองตัน เขตคลองเตย กรุงเทพมหานคร 10260
2533 Sukhumvit Road, Bangkok, Phrakhanong, Bangkok 10260
Phone: +66 2638 7000 Email: info.calibration@dksh.com Website: www.dksh.com/scientist-thailand

Delivering Growth - In Asia and Beyond.

เอกสารไม่ควบคุม
CAL-FM-C15-14: 06 Dec 2022

Reference standard equipment:

Equipment	Certificate no	Cal. date	Next Cal. date
Digital Thermometer with Probe	QR23-1073	2 May 23	2 May 24

Calibration Results:

Without Adjustment

Sensor Type: RTD

Channel: -

Diameter (mm) 4

Length (mm): 135

Immersion (mm): 110

Calibrate Point (°C)	STD. Reading (°C)	UUC. Reading (°C)	Correction of UUC (°C)	Uncertainty (± °C)
15.0	15.010	15.1	-0.090	0.076
25.0	25.006	25.1	-0.094	0.076
35.0	35.004	35.0	0.004	0.076

The End of Certificate

บริษัท ดีเคเอส อีเซีย จำกัด
DKSH Technology Limited
2533 ถนนสุขุมวิท แขวงคลองตัน เขตคลองเตย กรุงเทพมหานคร 10260
2533 Sukhumvit Road, Bangkok, Phrakhanong, Bangkok 10260
Phone: +66 2638 7000 Email: info.calibration@dksh.com Website: www.dksh.com/scientist-thailand

Delivering Growth - In Asia and Beyond.

เอกสารไม่ควบคุม
CAL-FM-C15-14: 06 Dec 2022

TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES & EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 16, SUKUMVIT, SUKUMVIT BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert.No.: 24MM292
Page.: 1 of 3

Equipment : Electronic Balance
Manufacturer : Mettler Toledo
Model : AB204-S/FACT
Serial No. : 1129361010
ID No. : UAE.WAS.002/2552
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangkok, Phrakhanong,
Bangkok 10260
Location : Balance Room (108)
Received order : 11 May 2024
Calibration Date : 11 May 2024
Ambient Temperature : 15 °C to 40 °C
Relative Humidity : 30 % to 90 %
Calibrated by : Khit Ruttanaprapachai
Approved by :
() Ponpan Paipim
() Suwit Imjai
(✓) Kunchit Promprat

Issue Date : 15 May 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3: Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : Electronic Balance
Condition As-Received : Used Item
Reference : 2405-01680C-1
Procedure used :-

Cert.No.: 24MM292
Page: 2 of 3

Calibration were conducted using in-house calibration procedure CP-OB01 based on UKAS LAB 14
according to direct measurement method against standard weight.

Condition of this result of calibration

1. Reference standard instruments:-

Instruments	Model	Serial No.	ID No.	Test report No.	Due date
1) Standard Weight Set (E2)	15884	24053	70RC007	MM-0013-24	25 Jan 2026

- This certificate is valid only to the item calibrated on date and place of calibration.
- This result of calibration was made on requested at the point specified by customer.
- This certificate is not certified for any commercial transaction.
- This certificate is traceable to the International System of Unit.

Result of calibration () Without Adjustment (*) After Adjustment by Internal Calibration

Range capacity : 0 g to 220 g Resolution 0.0001 g

Before Adjustment :

Applied Weight	Balance Reading	Correction	Measurement Uncertainty	Coverage Factor
(g)	(g)	(g)	(± mg)	(k)
100	100.0000	0.0000	0.19	2.03
200	200.0006	-0.0006	0.30	2

After Adjustment :

1. Determination of the standard deviation of weighing machine (n = 10)


Applied Weight	Standard Deviation of Reading (g)
(g)	
100	0.00007
200	0.00005

เอกสารไม่ควบคุม



Equipment : Electronic Balance
Condition As-Received : Used Item
Reference : 2405-0168OC-1
Result of calibration

Cert.No.: 24MM292
Page: 3 of 3

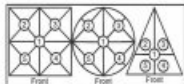

Maximum difference between off-center and central loading

Position 1 (g)	Position 2 (g)	Position 3 (g)	Position 4 (g)	Position 5 (g)	(g)
-0.0004	-0.0004	-0.0003	-0.0003	-0.0004	0.0001

3. Departure from nominal value

Applied Weight (g)	Balance Reading (g)	Correction (g)	Measurement Uncertainty (± mg)	Coverage Factor (k)
Unload	0.0000	0.0000	0.15	2.13
0.01	0.0100	0.0000	0.15	2.13
0.05	0.0500	0.0000	0.15	2.13
0.1	0.1000	0.0000	0.15	2.13
0.5	0.5000	0.0000	0.15	2.13
1	1.0000	0.0000	0.15	2.13
10	10.0000	0.0000	0.15	2.11
50	49.9999	+0.0001	0.17	2.06
100	99.9999	+0.0001	0.19	2.03
150	149.9998	+0.0002	0.29	2
200	199.9990	+0.0010	0.30	2

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k , providing a level of confidence of approximately 95 %.

-odo-

เอกสารไม่ควบคุม



Request No. 25-67 / 0275

MTC. ACL.No. 358 / 67

CALIBRATION CERTIFICATE

NOMENCLATURE : 1. Atomic Absorption Spectrophotometer "Agilent Technologies"

Model AA240FS, Serial No. MY13160001

2. Working standard solution "Inorganic Ventures"

Multi Analyte Custom Grade Solution, Lot No. S2-MEB675610

SUBMITTED BY : United Analyst and Engineering Consultant Co., Ltd.

3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Prakanong, Bangkok 10260

CALIBRATION PROCEDURE : 1. Performance Verification of Atomic Absorption Spectrophotometer
(WI-500-02-30)

2. Estimation Uncertainty of Measurement in Analytical Chemistry (QP-513)

CALIBRATION RANGE: 0.02, 0.10, 0.30, 0.50, 0.70 mg/L at 228.8 nm.Cd, 0.10, 0.20, 0.30, 0.50, 0.70 mg/L at 357.9 nm.Cr,
0.05, 0.10, 0.30, 0.50, 0.70 mg/L at 324.7 nm.Cu, 0.10, 0.30, 0.50, 0.70, 1.00 mg/L at 248.3 nm.Fe,
0.20, 0.50, 0.70, 1.00, 1.50 mg/L at 217.0 nm.Pb, 0.05, 0.10, 0.30, 0.50, 0.70 mg/L at 279.5 nm.Mn,
0.10, 0.30, 0.50, 0.70, 1.00 mg/L at 232.0 nm.Ni, 0.05, 0.10, 0.30, 0.50, 0.70 mg/L at 213.9 nm.Zn

CALIBRATION DATE : 2 February 2024

REFERENCE MATERIAL : Traceable to NIST "Agilent Technologies", "CARLO ERBA"

Cadmium Lot No. 0006589926, Chromium Lot No. 0112384886, Copper Batch No. T117098A, Iron Batch No. T126087A,
Lead Lot No. 1227873, Manganese Batch No. T109228A, Nickel Batch No. T270178A, Zinc Batch No. T820140A

AMBIENT CONDITIONS : Temperature 25 ± 5 °C. Relative humidity 50 ± 20 %

The Atomic Absorption Spectrophotometer has been calibrated against Reference Material traceable to National Institute of Standards and Technology (NIST) by The Analytical Chemistry Laboratory. The results are attached herewith.

Calibrated by Atipat

(Mr. Atipat Ratana)

Approved by S. S. S. S.

(Miss Sutadtha Deavong)

Director of Analytical Chemistry Laboratory

Ref. 2015267020100454001

Issued Date : 11 March 2024

The results relate only to the items tested/calibrated or value assigned.
Advertising the Report/Certificate and publicity of the results except in full are prohibited unless written permission is obtained from the governor of TISTR.

FMBL.MTC.002 Rev.4

Head Office
35 Mu 3 Tambon Khlong Ha, Amphoe Khlong Luang,
Changwat Pathumthani 12120, Thailand
Tel. (66) 0 2577 9000
Fax. (66) 0 2577 9009
E-mail : numpai@tistr.or.th Website:www.tistr.or.th

Office/Laboratory
Soi 1C, Bangpoo Industrial Estate, Sukhumvit Road,
Amphoe Muang, Changwat Samutprakan 10280, Thailand
Tel. (66) 0 2323 1672-80 ext. 115, 116
Fax. (66) 0 2323 9165
E-mail : mtg@tistr.or.th

Office
196 Phahonyothin Road, Chatuchak, Bangkok 10900,
Thailand
Tel. (66) 0 2577 9000
Fax. (66) 0 2577 9009
E-mail : sumalee@tistr.or.th

เอกสารไม่ควบคุม



Request No. 25-67 / 0275

1 / 5

MTC. ACL. No. 358 / 67

CALIBRATION DATA

1. Noise Level

Element	Cd	Cr	Cu	Fe	Pb	Mn	Ni	Zn
Absorbance	0.0006	0.0004	-0.0003	0.0001	-0.0011	-0.0005	0.0008	0.0004
	0.001	0.0017	-0.0009	0.0008	0.0001	0.0002	-0.0003	0.0007
	0.0006	0.0017	-0.0020	0.0005	0.0005	0.0004	0.0013	0.0014
	0.0001	0.0018	-0.0007	0.0005	0.0004	-0.0003	-0.0001	0.0010
	-0.0001	0.0019	-0.0014	0.0003	0.0010	0.0000	0.0002	-0.0001
	0.0011	0.0014	-0.0017	0.0009	-0.0008	0.0004	0.0006	0.0010
	-0.0002	0.0015	-0.0015	0.0003	0.0002	-0.0008	0.0009	0.0013
	0.0006	0.0012	-0.0001	0.0006	0.0008	0.0001	-0.0002	0.0013
	0.0008	0.0009	-0.0003	0.0003	0.0005	0.0002	0.0001	0.0007
	0.0012	0.0011	-0.0012	0.0008	0.0003	0.0004	0.0004	0.0013
	0.0003	0.0015	-0.0019	0.0001	-0.0002	0.0000	-0.0003	0.0003
	0.0005	0.0017	-0.0019	-0.0007	0.0000	-0.0007	0.0005	0.0005
	-0.0006	0.0016	0.0000	0.0006	-0.0001	0.0013	0.0006	0.0010
	0.0003	0.0011	-0.0002	0.0001	-0.0007	0.0009	0.0009	0.0002
	0.0003	0.0012	-0.0011	0.0007	-0.0003	0.0003	0.0010	0.0009
	0.0004	0.0018	-0.0016	-0.0004	-0.0006	0.0008	0.0007	0.0007
	-0.0001	0.0018	-0.0018	0.0013	-0.0006	-0.0001	0.0014	0.0006
	0.0003	0.0017	-0.0001	0.0001	-0.0012	-0.0004	0.0001	0.0002
	0.0010	0.0018	-0.0007	0.0003	-0.0005	-0.0002	0.001	0.0003
	0.0004	0.0019	-0.0008	-0.0001	-0.0004	0.0003	0.0002	0.0008
Average Absorbance	0.000	0.001	-0.001	0.000	0.000	0.000	0.000	0.001

Continue 2 / 5

INDUSTRIAL METROLOGY AND TESTING SERVICE CENTRE

The results relate only to the items tested/calibrated or value assigned.
Advertising the Report/Certificate and publicity of the results except in full are prohibited unless written permission is obtained from the governor of TISTR.

FMBL.MTC.002 Rev.4

Head Office
35 Mu 3 Tambon Khlong Ha, Amphoe Khlong Luang,
Changwat Pathumthani 12120, Thailand
Tel. (66) 0 2577 9000
Fax. (66) 0 2577 9009
E-mail : numpai@tistr.or.th Website:www.tistr.or.th

Office/Laboratory
Soi 1C, Bangpoo Industrial Estate, Sukhumvit Road,
Amphoe Muang, Changwat Samutprakan 10280, Thailand
Tel. (66) 0 2323 1672-80 ext. 115, 116
Fax. (66) 0 2323 9165
E-mail : mtg@tistr.or.th

Office
196 Phahonyothin Road, Chatuchak, Bangkok 10900,
Thailand
Tel. (66) 0 2577 9000
Fax. (66) 0 2577 9009
E-mail : sumalee@tistr.or.th

เอกสารไม่ควบคุม

Request No. 25-67 / 0275

2 / 5

MTC. ACL. No. 358 / 67

2. Precision

Element	Conc. (mg/L)	Absorbance										Ave. Abs.	SD	WRSD
Cd	0.02	0.0078	0.0076	0.0069	0.0075	0.0071	0.0070	0.0076	0.0074	0.0077	0.0067	0.007	0.0004	5.15
	0.30	0.1008	0.1007	0.0999	0.0997	0.1000	0.0996	0.1008	0.1002	0.1005	0.0999	0.100	0.0005	0.46
	0.70	0.2301	0.2306	0.2277	0.2305	0.2310	0.2295	0.2290	0.2293	0.2305	0.2296	0.230	0.0010	0.42
Cr	0.10	0.0094	0.0093	0.0093	0.0098	0.0094	0.0095	0.0090	0.0090	0.0094	0.0090	0.009	0.0003	2.75
	0.30	0.0241	0.0236	0.0221	0.0238	0.0231	0.0226	0.0231	0.0225	0.0230	0.0231	0.023	0.0006	2.75
	0.70	0.0500	0.0500	0.0500	0.0524	0.0499	0.0511	0.0509	0.0512	0.0515	0.0504	0.051	0.0008	1.63
Cu	0.05	0.0061	0.0062	0.0064	0.0061	0.0069	0.0061	0.0062	0.0064	0.0061	0.0061	0.006	0.0003	5.00
	0.30	0.0419	0.0411	0.0402	0.0407	0.0405	0.0404	0.0399	0.0400	0.0399	0.0400	0.040	0.0006	1.58
	0.70	0.0960	0.0960	0.0960	0.0959	0.0947	0.0955	0.0952	0.0951	0.0955	0.0956	0.095	0.0005	0.48
Fe	0.10	0.0096	0.0101	0.0103	0.0100	0.0099	0.0096	0.0106	0.0099	0.0105	0.0102	0.010	0.0003	3.38
	0.50	0.0424	0.0415	0.0428	0.0427	0.0421	0.0426	0.0413	0.0430	0.0421	0.0419	0.042	0.0006	1.33
	1.00	0.0830	0.0839	0.0847	0.0834	0.0832	0.0820	0.0839	0.0838	0.0837	0.0845	0.084	0.0008	0.92
Pb	0.20	0.0078	0.0074	0.0078	0.0078	0.0076	0.0078	0.0077	0.0078	0.0078	0.0077	0.008	0.0001	1.71
	0.70	0.0278	0.0273	0.0271	0.0267	0.0270	0.0264	0.0274	0.0273	0.0269	0.0269	0.027	0.0004	1.45
	1.50	0.0551	0.0548	0.0552	0.0555	0.0547	0.0546	0.0544	0.0549	0.0547	0.0547	0.055	0.0004	0.64
Mn	0.05	0.0116	0.0107	0.0110	0.0103	0.0108	0.0108	0.0112	0.0107	0.0109	0.0108	0.011	0.0003	3.15
	0.30	0.0650	0.0649	0.0649	0.0651	0.0646	0.0646	0.0649	0.0646	0.0640	0.0648	0.065	0.0003	0.48
	0.70	0.1463	0.1465	0.1459	0.1471	0.1475	0.1474	0.1487	0.1473	0.1462	0.1468	0.147	0.0008	0.56
Ni	0.10	0.0095	0.0100	0.0096	0.0103	0.0102	0.0096	0.0100	0.0095	0.0097	0.0096	0.010	0.0003	3.04
	0.50	0.0443	0.0433	0.0438	0.0444	0.0430	0.0437	0.0444	0.0437	0.0438	0.0434	0.044	0.0005	1.09
	1.00	0.0812	0.0820	0.0834	0.0829	0.0818	0.0829	0.0831	0.0835	0.0816	0.0819	0.082	0.0008	0.99
Zn	0.05	0.0374	0.0377	0.0373	0.0377	0.0374	0.0377	0.0373	0.0371	0.0374	0.0371	0.037	0.0002	0.61
	0.30	0.1985	0.1993	0.1975	0.1992	0.1979	0.1988	0.1995	0.1985	0.1974	0.2004	0.199	0.0009	0.47
	0.70	0.4027	0.4031	0.4019	0.4021	0.4023	0.3981	0.4042	0.4025	0.3993	0.3997	0.402	0.0019	0.48

Continue 3 / 5

INDUSTRIAL METROLOGY AND TESTING SERVICE CENTRE

The results relate only to the items tested/calibrated or value assigned.
Advertising the Report/Certificate and publicity of the results except in full are prohibited unless written permission is obtained from the governor of TISTR.

FMBL.MTC.002 Rev.4

Head Office
35 Mu 3 Tambon Khlong Ha, Amphoe Khlong Luang,
Changwat Pathumthani 12120, Thailand
Tel. (66) 0 2577 9000
Fax. (66) 0 2577 9009
E-mail : numpai@tistr.or.th Website:www.tistr.or.th

Office/Laboratory
Soi 1C, Bangpoo Industrial Estate, Sukhumvit Road,
Amphoe Muang, Changwat Samutprakan 10280, Thailand
Tel. (66) 0 2323 1672-80 ext. 115, 116
Fax. (66) 0 2323 9165
E-mail : mtg@tistr.or.th

Office
196 Phahonyothin Road, Chatuchak, Bangkok 10900,
Thailand
Tel. (66) 0 2577 9000
Fax. (66) 0 2577 9009
E-mail : sumalee@tistr.or.th

เอกสารไม่ควบคุม

3. Trueuess

3.1 Reading on wavelength- Cadmium(Cd) at 228.8 nm.

Element	Standard Value of RM (mg/l)	Reading (mg/l)	Error of Measurement (mg/l)	Error of Measurement (%)	Uncertainty (mg/l)
Cd	0.020	0.020	0.000	1.10	± 0.005
	0.301	0.301	0.000	0.11	± 0.005
	0.707	0.693	-0.013	1.85	± 0.008

3.2 Reading on wavelength- Chromium (Cr) at 357.9 nm.

Element	Standard Value of RM (mg/l)	Reading (mg/l)	Error of Measurement (mg/l)	Error of Measurement (%)	Uncertainty (mg/l)
Cr	0.1007	0.104	0.004	3.49	± 0.009
	0.3035	0.297	-0.006	2.11	± 0.012
	0.7071	0.685	-0.023	3.19	± 0.023

3.3 Reading on wavelength- Copper (Cu) at 324.7 nm.

Element	Standard Value of RM (mg/l)	Reading (mg/l)	Error of Measurement (mg/l)	Error of Measurement (%)	Uncertainty (mg/l)
Cu	0.051	0.047	-0.004	7.58	± 0.003
	0.303	0.296	-0.007	2.19	± 0.009
	0.704	0.698	-0.005	0.74	± 0.020

Continue 4 / 5

INDUSTRIAL METROLOGY AND TESTING SERVICE CENTRE

The results relate only to the items tested/calibrated or value assigned.

Advertising the Report/Certificate and publicity of the results except in full are prohibited unless written permission is obtained from the governor of TISTR.

FM/BLMTC.002 Rev.4

Head Office
35 Mu 3 Tambon Khlong Ha, Amphoe Khlong Luang,
Changwat Pathumthani 12120, Thailand
Tel. (66) 0 2577 9000
Fax. (66) 0 2577 9009
E-mail : numpai@tistr.or.th Website: www.tistr.or.th

Office/Laboratory
Sri 1C, Bangpoo Industrial Estate, Sukhumvit Road,
Amphoe Muang, Changwat Samutprakan 10280, Thailand
Tel. (66) 0 2323 1672-80 ext. 115, 116
Fax. (66) 0 2323 9165
E-mail : mtc@tistr.or.th

Office
196 Phahonyothin Road, Chatuchak, Bangkok 10900,
Thailand
Tel. (66) 0 2577 9000
Fax. (66) 0 2577 9009
E-mail : sumalee@tistr.or.th

เอกสารไม่ควบคุม

3.7 Reading on wavelength- Nickel (Ni) at 232.0 nm.

Element	Standard Value of RM (mg/l)	Reading (mg/l)	Error of Measurement (mg/l)	Error of Measurement (%)	Uncertainty (mg/l)
Ni	0.101	0.098	-0.003	2.90	± 0.013
	0.508	0.502	-0.006	1.16	± 0.018
	1.012	0.962	-0.051	5.02	± 0.032

3.8 Reading on wavelength- Zinc (Zn) at 213.9 nm.

Element	Standard Value of RM (mg/l)	Reading (mg/l)	Error of Measurement (mg/l)	Error of Measurement (%)	Uncertainty (mg/l)
Zn	0.050	0.045	-0.005	9.39	± 0.013
	0.303	0.324	0.021	7.04	± 0.013
	0.707	0.675	-0.032	4.52	± 0.019

Remark : The reported uncertainty is an expanded uncertainty calculated using a coverage factor of 2 (k = 2)
which gives a level of confidence of approximately 95%

Calibrated by Atipat
(Mr. Atipat Ratana)

Approved by Suladda
(Miss Suladda Deawtong)
Director of Analytical Chemistry Laboratory
Issued Date : 11 March 2024

INDUSTRIAL METROLOGY AND TESTING SERVICE CENTRE

End of Certificate

The results relate only to the items tested/calibrated or value assigned.

Advertising the Report/Certificate and publicity of the results except in full are prohibited unless written permission is obtained from the governor of TISTR.

FM/BLMTC.002 Rev.4

Head Office
35 Mu 3 Tambon Khlong Ha, Amphoe Khlong Luang,
Changwat Pathumthani 12120, Thailand
Tel. (66) 0 2577 9000
Fax. (66) 0 2577 9009
E-mail : numpai@tistr.or.th Website: www.tistr.or.th

Office/Laboratory
Sri 1C, Bangpoo Industrial Estate, Sukhumvit Road,
Amphoe Muang, Changwat Samutprakan 10280, Thailand
Tel. (66) 0 2323 1672-80 ext. 115, 116
Fax. (66) 0 2323 9165
E-mail : mtc@tistr.or.th

Office
196 Phahonyothin Road, Chatuchak, Bangkok 10900,
Thailand
Tel. (66) 0 2577 9000
Fax. (66) 0 2577 9009
E-mail : sumalee@tistr.or.th

เอกสารไม่ควบคุม

3.4 Reading on wavelength- Iron (Fe) at 248.3 nm.

Element	Standard Value of RM (mg/l)	Reading (mg/l)	Error of Measurement (mg/l)	Error of Measurement (%)	Uncertainty (mg/l)
Fe	0.100	0.104	0.005	4.60	± 0.014
	0.500	0.482	-0.018	3.55	± 0.016
	1.006	0.968	-0.038	3.75	± 0.029

3.5 Reading on wavelength- Lead (Pb) at 217.0 nm.

Element	Standard Value of RM (mg/l)	Reading (mg/l)	Error of Measurement (mg/l)	Error of Measurement (%)	Uncertainty (mg/l)
Pb	0.201	0.202	0.001	0.34	± 0.014
	0.706	0.719	0.012	1.73	± 0.030
	1.513	1.459	-0.054	3.57	± 0.061

3.6 Reading on wavelength- Manganese (Mn) at 279.5 nm.

Element	Standard Value of RM (mg/l)	Reading (mg/l)	Error of Measurement (mg/l)	Error of Measurement (%)	Uncertainty (mg/l)
Mn	0.0505	0.050	0.000	0.83	± 0.005
	0.3031	0.306	0.003	1.12	± 0.007
	0.7023	0.698	-0.004	0.62	± 0.014

Continue 5 / 5

INDUSTRIAL METROLOGY AND TESTING SERVICE CENTRE

The results relate only to the items tested/calibrated or value assigned.

Advertising the Report/Certificate and publicity of the results except in full are prohibited unless written permission is obtained from the governor of TISTR.

FM/BLMTC.002 Rev.4

Head Office
35 Mu 3 Tambon Khlong Ha, Amphoe Khlong Luang,
Changwat Pathumthani 12120, Thailand
Tel. (66) 0 2577 9000
Fax. (66) 0 2577 9009
E-mail : numpai@tistr.or.th Website: www.tistr.or.th

Office/Laboratory
Sri 1C, Bangpoo Industrial Estate, Sukhumvit Road,
Amphoe Muang, Changwat Samutprakan 10280, Thailand
Tel. (66) 0 2323 1672-80 ext. 115, 116
Fax. (66) 0 2323 9165
E-mail : mtc@tistr.or.th

Office
196 Phahonyothin Road, Chatuchak, Bangkok 10900,
Thailand
Tel. (66) 0 2577 9000
Fax. (66) 0 2577 9009
E-mail : sumalee@tistr.or.th

เอกสารไม่ควบคุม

Agilent CrossLab Start Up Services

Agilent 5100 5110 ICP-OES
Preventive Maintenance

Agilent Preventive Maintenance provides factory recommended service for your analytical instruments to assure reliable operation and the accuracy of your results

Delivered by highly trained and certified service engineers using genuine Agilent parts and supplies, Agilent Preventive Maintenance provides what you need to reduce unplanned downtime and keep your systems operating at their peak performance.

This checklist is used as a guide for completing the preventive maintenance tasks. A signed copy of this checklist is provided for your records.

Introduction

Customer Information

- Customers should provide all necessary operating supplies upon request of the engineer.
- A customer representative should be available to the engineer while performing the preventive maintenance procedures. Customers are responsible for regular maintenance and are encouraged to observe the service representative.
- Any parts not included in the Parts Lists section of this document are not part of the recommended Preventive Maintenance service nor are they included in the price of this service.
- If a system requires the use of extra or special procedures and/or parts for the maintenance service, then these must be ordered separately and charged as a repair, which may incur additional costs.
- For customers using HF applications, the instrument should be returned to its standard sample introduction system.

Important Customer Web Links

- To access **Agilent University**, visit <http://www.agilent.com/crosslab/university/> to learn about training options, which include online, classroom and onsite delivery. A training specialist can work directly with you to help determine your best options.
- To access the **Agilent Resource Center** web page, visit <https://www.agilent.com/en-us/agilentresources>. The following information topics are available:
 - Sample Prep and Containment
 - Chemical Standards
 - Analysis
 - Service and Support
 - Application Workflows
- The **Agilent Community** is an excellent place to get answers, collaborate with others about applications and Agilent products, and find in-depth documents and videos relevant to Agilent technologies. Visit <https://community.agilent.com/welcome>
- Videos about specific preparation requirements for your instrument can be found by searching the **Agilent YouTube** channel at <https://www.youtube.com/user/agilent>
- Need to place a service call?** [Flexible Repair Options | Agilent](#)

Service Engineer's Responsibilities

- Contact the customer and ensure that all necessary supplies are available before the preventive maintenance visit.
- Only select those pages that relate to the system or module being serviced.
- Complete empty fields with the relevant information.
- Complete the relevant checkboxes in the checklist using either a "X" or tick mark "✓".
- Check **"Service not applicable"** check boxes to indicate services/tasks not delivered, as appropriate.
- Complete the Preventive Maintenance services in the most logical order relevant to the individual system service in the order of the tasks listed.
- Complete the **Service Review** section together with the customer.
- Complete the fields for page numbers at the foot of each selected page.
- Add relevant page numbers to selected pages and complete the total number of pages field in the Service Completion section.
- Ask the customer to sign the Service Verification section including the customer's and your signature.

Instrument Maintenance

System Information

- ☐ Check this box if an instrument configuration report is attached instead of completing the table.

Instrument System Name and ID	5110 VDV ICP-OES
Instrument System Site and Location	United Analytical and Engineering Consultant

List System Component Product Numbers	List the Serial Numbers of each Component
1. G 3000	77 16030001
2.	
3.	
4.	
5.	
6.	
7.	
8.	
9.	

ICP-OES Configuration Table	Circle the type or write in the type if other
Nebulizer Type	SeaSpray (One Neb) Conical Other
Spray Chamber	Cyclonic Single Pass (Cyclonic Double Pass) Other
Torch	Radial (Dual View) Other
Torch Type	One Piece (Semi Demountable) Fully Demountable Other
Injector Diameter	2.4mm (1.8mm) 1.4mm 0.8mm Other
Injector Material	Quartz Ceramic Other

Preparation

- ☒ Discuss any specific issues with the customer before starting.
- ☒ Review the instrument logbook for recorded problems and comments.
- ☒ Save instrument control settings before starting the procedure.
- ☒ Perform a general inspection of the system for cleanliness.
- ☒ Check for proper installation of parts, assemblies, sensors etc.
- ☒ Check system for required installation of components and implementation of Service Notes
- ☒ Check for required firmware/software updates and verify with customers if they would like them installed.
- ☐ For HF application systems, if standard sample introduction system was not installed, ask the customer to install it. *vi h*
- ☒ Ask the customer to remove any samples from the ICP-OES sample introduction area, auto sampler or around the ICP-OES.

Preventive Maintenance Procedures

Record Pre-PM instrument performance

- ☒ Run Instrument Performance test.
- ☒ Record results in Instrument Performance Test Results Table – Pre-PM.

Clean and inspect ICP-OES system

- ☒ Look for any obvious external damage or problems.
- ☒ Inspect water cooling hoses, gas lines and power cord for excessive wear or damage.
- ☒ Perform a general internal inspection of the system for excessive dust accumulation, clean if necessary.
- ☒ Inspect sample introduction components and record any required maintenance in the Service Engineer Comments and notify the customer as the required actions required.
- ☒ Record the instrument operating conditions in the ICP-OES Status Results Table.
- ☒ Replace the polychromator purge filter.
- ☒ Replace the radial pre-optics window
- ☒ Replace the axial pre-optics window for SVDV and VDV instruments.
- ☒ Check exhaust flow for the correct positive extraction at the exhaust duct to insure they meet minimum specifications.
- ☒ Replace air inlet dust filter.
- ☐ Replace high capacity air inlet dust filter element if installed. *vi h*
- ☒ Remove and clean instrument water inlet filter.

Agilent Water Recirculator

- ☐ **Service not applicable**
- ☒ Drain cooling fluid and remove any particles from the chiller reservoir
- ☒ Remove, clean and reinstall water inlet metal mesh filter if present.
- ☒ Re fill with Agilent Cool Clear cooling fluid.
- ☒ Clean the cooling system Air filter and the condenser.

SPS 3 Auto Sampler

- ☒ **Service not applicable**
- ☐ Power cycle the autosampler and verify successful initialization.
- ☐ Inspect X and Z axis belts for wear. Replace is necessary.
- ☐ Clean X and Z axis slide shafts.
- ☐ Using customer's racks and the Agilent software move the sample probe to the 4 outermost corners and rinse port, ensure that the probe is approximately centered in the vial.

SPS 4 Auto sampler

- ☒ **Service not applicable**
- ☐ Clean the spill tray, rack location mat, end frames and chassis with a damp soft cloth and diluted mild detergent.
- ☐ Clean the auto sampler cover panels, if cover kit is installed, with domestic window cleaner.
- ☐ Check the X-axis and Z-axis drive belts for cracks, splits, damaged teeth, excessive fraying, color changes or degradation from fumes.
- ☐ Check the X-axis, Theta-axis and Z-axis FFC cables for cracks, incorrect positioning, damaged edges or damaged connectors.
- ☐ Pump Tubing Replacement. Replace peristaltic pump tubing. Replace all tubing that goes from the rinse station to the pump and from the pump to the waste/rinse bottles
- ☐ Test using customer's tray and move the sample probe to the sample vial 1, wash vial and rinse port and ensure that the probe is centered in the vial. If not use calibration wizard and calibrate the position.

AVS 4, 6, 7 Advanced Valve System

- ☒ **Service not applicable**
- ☐ Replace valve rotor seal
- ☐ Check fittings for signs of leaks
- ☐ Check tubing including autosampler tubing for kinks or excessive wear
- ☐ Check high flow pump for signs of leaks

ICP-OES adjustment

- ☒ Check position of Zn peak, adjust if required.
- ☒ Check Argon Ratio, adjust to specified value if required.
- ☒ Perform Detector Calibration.
- ☒ Perform Instrument Calibration.

Record Post-PM instrument performance

- ☒ Run Instrument Performance test.
- ☒ Record results in Instrument Performance Test Results Table - Post PM.
- ☒ For systems using ICP Expert version 7.3 and above, run the following instrument tests
 - ☒ Subsystem Communications Test
 - ☒ Air Flow
 - ☒ Water Flow
 - ☒ Gas Flows
 - ☒ RF Generator
 - ☒ Camera Test
 - ☒ Optics Test
 - ☒ Nebulizer Test
- ☒ Record the result in the Instrument Test Results Table

Restore Instrument

- ☐ For HF applications, ask the customer to reinstall their sample introduction system. ^{N/A}
- ☒ Leave system in an idle state: on and purging.
- ☒ Guidance: If the PM service is performed prior to a qualification service, then use the qualification procedure as a guide for final instrument set up and checkout.

Service Review

- ☒ Attach available reports/printouts of all tests to this documentation.
- ☒ Record the Preventive Maintenance service activity in the customer's records/logbook.
- ☒ Record the PM event in the Smart Alerts logbook, if applicable.
- ☒ Update/reset instrument maintenance counters as appropriate.
- ☒ Affix the PM sticker to the system or instrument logbook based on the customer's request.
- ☒ Complete the Service Engineer Comments section if there are additional comments.
- ☒ Review this service, parts replaced, and test results obtained with the customer.
- ☒ If the instrument firmware was updated, record the details of the change in the Service Engineer's Comments box. Systems in a compliant environment may need additional documentation.
- ☒ Complete the Signature Page with both Service Engineer and Customer signatures.

Test Results

Instrument Performance Test Results Table

Note: These measurements do not form part of any specification and are for reference only.

	Pre PM Sensitivity Check		Post PM Sensitivity Check	
	Radial	Axial *	Radial	Axial*
Zn 213.857 nm SRBR	1500. %	2213.4	4124. %	6965.9
Mn 257.610 nm SRBR	3915.0	7492.1	13017. %	31121.6
Al 395.152 nm SBR	9.9	10.7	9.9	21.1
K 766.491 nm SBR	5.7	23.1	4. %	45.3

* Axial result is not applicable for G8016AA, G8012AA Radial View instruments.

Instrument Test Results Table

Note: The Instrument Test results are for systems using ICP Expert version 7.3 and above only.

Instrument Test	Result
Subsystem Communications Test	Pass
Air Flow	Pass
Water Flow	Pass
Gas Flows	Pass
RF Generator	Pass
Camera Test	Pass
Optics Test	Pass
Nebulizer test	Pass

เอกสารไม่ควบคุม

เอกสารไม่ควบคุม

ICP-OES Status Results Table

Note: These measurements do not form part of any specification and are for reference only.

Measurement	Standby Mode	Plasma On
Mains Voltage	231.411 VAC	226.931 VAC
Mains Current	0.0 % A	0.105 A
Instrument Temperature	22.1 °C	23.5 °C
RF Air Flow (sensor speed)	14.0 Hz	19.0 Hz
Plasma Exhaust Temperature	No measurement	63. % °C
Water Flow Oscillator	No measurement	1.34 L/min
Water Flow Detector	0.86 L/min	0.81 L/min
Water Inlet Temperature	19.7 °C	19.9 °C
Polychromator Temperature	35.0 °C	35.0 °C
CCD Temperature	-40.1 °C	-39.6 °C
Thermal Stabilizer	35.0 °C	35.0 °C
Argon Supply Pressure	649.92 kPa	591.55 kPa
Purge Gas Supply Pressure*1	646.66 kPa	612.41 kPa
Option Gas Supply Pressure*1	- kPa	- kPa
Nebulizer Flow	No measurement	0.90 L/min
Nebulizer Back Pressure	No measurement	158.43 kPa
Plasma Gas Flow	No measurement	11.91 L/min
Auxiliary Gas Flow	No measurement	1.00 L/min
RF Power	No measurement	1204.7 W
RF Supply Current	No measurement	7.658 A
RF Supply Voltage	No measurement	204.417 V

*1 If option installed

Consumed PM Parts

Part Description	Part Number	Product or Model# where used	Quantity consumed
Axial Pre-Optic Window	G8010-68014	G8010A, G8011A, G8014A/G8015A	1
Radial Pre-Optic Window	G8010-68015	All	1
Agilent Cool Clear Coolant Fluid	5799-0037	Agilent Water Recirculator	-
Purge Gas Filter	G8010-60136	All	1
Air inlet filter	G8000-68002	All	1
High Capacity Air Filter	G8010-60189	Optional	-
Rotor seal for 6-7 port valve for AVS6/7	G8494-60002	G8494A/G8495	-
Rotor seal for 4 port valve for AVS4	G8493-60002	G8493A	-
Rinse solution to rinse station 2.5mm id x 1m	G8410-80123	SPS 4	-
Barb connector 2.5mm-1.5mm ID	G8410-80124	SPS 4	-
PVC waste tubing, 8mm od x 5mm id, 2m	G8410-80122	SPS 4	-
Additional Parts may be required from engineer's stock:			
X axis drive belt	5410047500	SPS 3	-
Z axis drive belt	5410047400	SPS 3	-
Peristaltic pump tubing, PVC SolveFlex, 3 bridged,	3710049000	SPS 4	-

Consumed Parts Reference (Purchased by customer, not included as part of PM)

☒ Section Not Applicable

Part Description	Part Number	Product or Model# where used	Quantity consumed

เอกสารไม่ควบคุม

เอกสารไม่ควบคุม

Signature Page

Service Engineer Comments (optional)

If there are any specific points you wish to note as part of performing the installation or other items of interest for the customer, please write in this box.

Service Verification

Service Request Number:

6003193100

Date Service Completed:

04 Nov 2024

Service Engineer Name:

Kanyakorn S.

Customer Name:

Aphorn Onkong

Service Engineer Signature:

Kanyakorn S.

Customer Signature:

Aphorn Onkong

Total number of pages in this document:

14

เอกสารไม่ควบคุม

Report Summary

Instrument Model Agilent 5100/5110 VDV ICP-OES
Instrument ID G8011A/G8015A
Instrument Serial Number MY18030001
Software Version 7.3.1.9507
Firmware Version 3442
Tested By Pre Test_PM_Kanyakorn S.
Test Completed On 11/4/2024 9:19:10 AM

Result Summary

Subsystem Communications Test	Skipped
Air Flow Test	Skipped
Water Flow Test	Skipped
Gas Flows Test	Skipped
RF Generator Test	Skipped
Camera Test	Skipped
Optics Test	Skipped
Advanced Valve System Test	Skipped
Resolution Test	Pass
Sensitivity Test	Fail
Precision Test	Pass

เอกสารไม่ควบคุม

Resolution Test

Pass

Element Wavelength	Specification	Width
N (174.213 nm)	≤ 9.40	6.98
As (188.980 nm)	≤ 8.20	6.17
C (193.027 nm)	≤ 11.50	8.30
Mo (202.032 nm)	≤ 8.20	6.38
Cr (206.158 nm)	≤ 13.40	8.88
Zn (213.857 nm)	≤ 8.70	6.80
Pb (220.353 nm)	≤ 9.50	7.09
Co (228.615 nm)	≤ 17.20	11.67
Ba (230.424 nm)	≤ 9.40	7.20
Mn (257.610 nm)	≤ 13.30	9.43
Mn (259.588 nm)	≤ 20.30	14.11
Cr (267.716 nm)	≤ 11.00	8.04
Cu (324.754 nm)	≤ 25.00	18.97
Cu (327.395 nm)	≤ 14.20	11.23
Sr (338.071 nm)	≤ 33.50	24.30
Ba (455.403 nm)	≤ 44.00	33.47
Sr (460.733 nm)	≤ 36.00	17.23
Ba (493.408 nm)	≤ 35.00	25.37
Ba (614.171 nm)	≤ 42.00	25.54
Ar (675.283 nm)	≤ 74.00	56.51
K (766.491 nm)	≤ 80.00	65.95

เอกสารไม่ควบคุม

Sensitivity Test

Fail

Radial					
Element Wavelength	Specification	Method	Ratio	Standard	Blank
As (188.980 nm)	≥ 46.0	SRBR	104.1	793.0	50.8
Se (196.026 nm)	≥ 41.0	SRBR	87.6	882.0	79.7
Zn (213.857 nm)	≥ 1421.0	SRBR	1500.8	41823.3	749.0
Pb (220.353 nm)	≥ 46.0	SRBR	170.7	2432.0	174.9
Mn (257.610 nm)	≥ 3518.0	SRBR	3915.0	264700.2	4420.0
Al (396.152 nm)	≥ 3.4	SBR	7.7	48454.6	5563.2
Ba (493.408 nm)	≥ 34.0	SBR	45.9	1966719.7	41903.8
K (766.491 nm)	≥ 1.8	SBR	5.7	99038.2	14687.7
Axial					
Element Wavelength	Specification	Method	Ratio	Standard	Blank
As (188.980 nm)	≥ 208.0	SRBR	126.5	1498.8	119.0
Se (196.026 nm)	≥ 159.0	SRBR	112.0	1773.6	157.8
Zn (213.857 nm)	≥ 234.0	SRBR	466.0	6784.2	199.7
Cd (214.439 nm)	≥ 4227.0	SRBR	1919.3	68724.6	1236.4
Pb (220.353 nm)	≥ 320.0	SRBR	332.5	7929.5	499.0
Mn (257.610 nm)	≥ 10625.0	SRBR	7492.2	991238.3	16911.7
Cr (267.716 nm)	≥ 1048.0	SRBR	2254.6	129706.6	3150.9
Cu (324.754 nm)	≥ 19.0	SBR	26.9	290748.3	10407.5
Al (396.152 nm)	≥ 6.0	SBR	10.7	211329.2	18005.0
Ba (493.408 nm)	≥ 60.0	SBR	49.3	6956460.4	138335.9
K (766.491 nm)	≥ 24.0	SBR	26.1	1395190.2	47996.2

เอกสารไม่ควบคุม

Precision Test			Pass
Radial			
Element Wavelength	Specification	Measured Value % RSD	
As (188.980 nm)	≤ 2.60	0.73	
Se (196.026 nm)	≤ 2.60	0.95	
Zn (213.857 nm)	≤ 1.50	0.31	
Pb (220.353 nm)	≤ 2.60	0.73	
Mn (257.610 nm)	≤ 1.50	0.39	
Al (396.152 nm)	≤ 1.50	0.39	
Ba (493.408 nm)	≤ 1.50	0.87	
K (766.491 nm)	≤ 1.50	0.32	
Axial			
Element Wavelength	Specification	Measured Value % RSD	
As (188.980 nm)	≤ 1.50	1.21	
Se (196.026 nm)	≤ 1.50	0.84	
Zn (213.857 nm)	≤ 1.50	0.56	
Zn (213.857 nm)	≤ 1.50	0.98	
Cd (214.439 nm)	≤ 1.50	0.28	
Pb (220.353 nm)	≤ 1.50	0.51	
Mn (257.610 nm)	≤ 1.50	0.97	
Cr (267.716 nm)	≤ 1.50	0.22	
Cu (324.754 nm)	≤ 1.50	0.24	
Al (396.152 nm)	≤ 1.50	0.33	
Ba (493.408 nm)	≤ 1.50	0.40	
K (766.491 nm)	≤ 1.50	0.65	

Page 4 of 4

เอกสารไม่ควบคุม

Report Summary	
Instrument Model	Agilent 5100/5110 VDV ICP-OES
Instrument ID	G8011A/G8015A
Instrument Serial Number	MY18030001
Software Version	7.3.1.9507
Plate/Well	3442
Tested By	Post Test_PM_Kanyakorn S.
Test Completed On	11/4/2024 11:07:24 AM
Result Summary	
Subsystem Communications Test	Pass
Air Flow Test	Skipped
Water Flow Test	Skipped
Gas Flows Test	Skipped
Low Pressure Vapour	Skipped
RF Generator Test	Skipped
Camera Test	Skipped
Optics Test	Pass
Advanced Valve System Test	Skipped
Resolution Test	Pass
Sensitivity Test	Fail
Precision Test	Pass
Subsystem Communications Test	Pass
Optics Test	
	Pass
Intensity	Radial 3184054 Axial 3177175
Wavelength	Radial 737.212 Axial 737.212

Page 1 of 4

เอกสารไม่ควบคุม

Resolution Test			Pass
Element Wavelength	Specification	Width	
N (174.213 nm)	≤ 9.40	6.97	
As (188.980 nm)	≤ 8.20	6.14	
C (193.027 nm)	≤ 11.50	8.33	
Mo (202.032 nm)	≤ 8.20	6.33	
Cr (267.138 nm)	≤ 13.40	9.06	
Zn (213.857 nm)	≤ 8.70	6.70	
Pb (220.353 nm)	≤ 9.50	7.03	
Co (228.615 nm)	≤ 17.20	11.72	
Ba (230.424 nm)	≤ 9.40	7.32	
Mn (257.610 nm)	≤ 13.30	9.44	
Mn (260.568 nm)	≤ 20.30	14.21	
Cr (267.716 nm)	≤ 11.00	7.94	
Cu (324.754 nm)	≤ 25.00	18.99	
Cu (327.395 nm)	≤ 14.20	11.27	
Sr (338.071 nm)	≤ 33.50	24.40	
Ba (455.403 nm)	≤ 44.00	33.50	
Sr (460.753 nm)	≤ 36.00	17.31	
Ba (493.408 nm)	≤ 36.00	25.44	
Ba (614.171 nm)	≤ 42.00	25.16	
Ar (675.283 nm)	≤ 74.00	56.15	
K (766.491 nm)	≤ 80.00	65.56	

Page 2 of 4

เอกสารไม่ควบคุม

Sensitivity Test						Fail
Radial						
Element Wavelength	Specification	Method	Ratio	Standard	Blank	
As (188.980 nm)	≥ 46.0	SRBR	130.6	977.1	50.4	
Se (196.026 nm)	≥ 41.0	SRBR	108.0	958.7	70.2	
Zn (213.857 nm)	≥ 1421.0	SRBR	4124.8	44037.7	113.4	
Pb (220.353 nm)	≥ 46.0	SRBR	207.2	2554.7	136.2	
Mn (257.610 nm)	≥ 3518.0	SRBR	13017.8	271846.6	434.7	
Al (396.152 nm)	≥ 3.4	SBR	9.7	50615.5	4717.0	
Ba (493.408 nm)	≥ 34.0	SBR	133.7	2069203.0	15359.3	
K (766.491 nm)	≥ 1.6	SBR	4.8	100199.5	17235.5	
Axial						
Element Wavelength	Specification	Method	Ratio	Standard	Blank	
As (188.980 nm)	≥ 208.0	SRBR	174.9	1566.7	73.0	
Se (196.026 nm)	≥ 159.0	SRBR	187.0	1863.4	110.2	
Zn (213.857 nm)	≥ 234.0	SRBR	740.9	6836.0	83.1	
Zn (213.857 nm)	≥ 1743.0	SRBR	5965.9	101568.1	211.7	
Cd (214.439 nm)	≥ 4227.0	SRBR	5781.0	72852.9	158.1	
Pb (220.353 nm)	≥ 320.0	SRBR	501.0	8464.3	267.7	
Mn (257.610 nm)	≥ 10625.0	SRBR	31121.6	1006637.8	1044.0	
Cr (267.716 nm)	≥ 1048.0	SRBR	4424.8	132202.9	880.8	
Cu (324.754 nm)	≥ 19.0	SBR	68.7	302907.8	4345.6	
Al (396.152 nm)	≥ 6.0	SBR	21.1	218771.0	9892.3	
Ba (493.408 nm)	≥ 60.0	SBR	250.6	7137380.9	26367.3	
K (766.491 nm)	≥ 24.0	SBR	45.3	1435050.6	31025.0	

Page 3 of 4

เอกสารไม่ควบคุม

Precision Test			Pass
Radial			
Element Wavelength	Specification	Measured Value % RSD	
As (188.980 nm)	≤ 2.60	0.81	
Se (196.026 nm)	≤ 2.60	0.98	
Zn (213.857 nm)	≤ 1.50	0.22	
Pb (220.353 nm)	≤ 2.60	0.37	
Mn (257.610 nm)	≤ 1.50	0.27	
Al (396.152 nm)	≤ 1.50	0.25	
Ba (493.408 nm)	≤ 1.50	0.53	
K (766.491 nm)	≤ 1.50	0.15	
Axial			
Element Wavelength	Specification	Measured Value % RSD	
As (188.980 nm)	≤ 1.50	0.81	
Se (196.026 nm)	≤ 1.50	0.65	
Zn (213.857 nm)	≤ 1.50	0.79	
Zn (213.857 nm)	≤ 1.50	0.81	
Cd (214.439 nm)	≤ 1.50	0.35	
Pb (220.353 nm)	≤ 1.50	0.33	
Mn (257.610 nm)	≤ 1.50	1.02	
Cr (267.716 nm)	≤ 1.50	0.32	
Cu (324.754 nm)	≤ 1.50	0.51	
Al (396.152 nm)	≤ 1.50	0.37	
Ba (493.408 nm)	≤ 1.50	0.68	
K (766.491 nm)	≤ 1.50	0.74	

Page 4 of 4

เอกสารไม่ควบคุม

Report Summary	
Instrument Model	Agilent 5100/5110 VDV ICP-OES
Instrument ID	G8011A/G8015A
Instrument Serial Number	MY18030001
Software Version	7.3.1.9507
Firmware Version	3442
Tested By	Post Test_PM_Kanyakorn S.
Test Completed On	11/4/2024 11:30:15 AM
Result Summary	
Subsystem Communications Test	Pass
Air Flow Test	Pass
Water Flow Test	Pass
Gas Flows Test	Pass
RF Generator Test	Pass
Camera Test	Pass
Optics Test	Skipped
Advanced Valve System Test	Skipped
Resolution Test	Skipped
Sensitivity Test	Skipped
Precision Test	Skipped
Subsystem Communications Test	Pass
Air Flow Test	Pass
30% Air Flow (relative speed)	15.00
75% Air Flow (relative speed)	19.00
Water Flow Test	Pass
RF Water Flow(L/min)	1.30
Camera Water Flow (L/min)	0.81
Water Inlet Temperature (°C)	20.55

Page 1 of 2

เอกสารไม่ควบคุม

Gas Flows Test			Pass		
Nebulizer Target Flow	Actual Flow	Back Pressure	Auxiliary Target Flow	Actual Flow	Back Pressure
0.70	0.70	154.65	2.00	2.00	110.92
Makeup Target Flow	Actual Flow	Back Pressure	Plasma Target Flow	Actual Flow	Back Pressure
2.00	2.00	115.38	16.00	17.97	21.48
RF Generator Test			Pass		
RF Power Supply Test	Passed				
RF Power Supply (V)	128.554				
RF Oscillator Test	Passed				
RF Oscillator Frequency (MHz)	25.834				
Work Coil Current (A)	44.660				
RF Power Supply Current (A)	1.999				
Camera Test			Pass		
	Integration Time (ms)	Standard Deviation	Status		
Electronic Offset Test	1000	5.228	Passed		
Dark Current Test	6000	1.168	Passed		
Array Test	5	0.024	Passed		
Linearity Test		0.118	Passed		

Page 2 of 2

เอกสารไม่ควบคุม

Report Summary	
Instrument Model	Agilent 5100/5110 VDV ICP-OES
Instrument ID	G8011A/G8015A
Instrument Serial Number	MY18030001
Software Version	7.3.1.9507
Firmware Version	3442
Tested By	change mirror
Test Completed On	11/5/2024 10:35:26 AM
Result Summary	
Subsystem Communications Test	Skipped
Air Flow Test	Skipped
Water Flow Test	Skipped
Gas Flows Test	Skipped
RF Generator Test	Skipped
Camera Test	Skipped
Optics Test	Skipped
Advanced Valve System Test	Skipped
Resolution Test	Pass
Sensitivity Test	Pass
Precision Test	Pass

Page 1 of 4

เอกสารไม่ควบคุม

Resolution Test			Pass
Element Wavelength	Specification	Width	
N (174.213 nm)	≤ 9.40	8.79	
As (188.980 nm)	≤ 8.20	5.80	
C (193.027 nm)	≤ 11.50	8.15	
Mo (202.032 nm)	≤ 8.20	5.90	
Cr (206.158 nm)	≤ 13.40	8.85	
Zn (213.857 nm)	≤ 8.70	6.77	
Pb (220.353 nm)	≤ 9.50	6.61	
Co (228.615 nm)	≤ 17.20	11.79	
Ba (230.424 nm)	≤ 9.40	7.25	
Mn (257.610 nm)	≤ 13.30	9.47	
Mn (260.568 nm)	≤ 20.30	14.50	
Cr (267.716 nm)	≤ 11.00	7.91	
Cu (324.754 nm)	≤ 25.00	18.72	
Cu (327.395 nm)	≤ 14.20	11.09	
Sr (338.071 nm)	≤ 33.50	25.39	
Ba (455.403 nm)	≤ 44.00	33.09	
Sr (460.793 nm)	≤ 36.00	18.54	
Ba (493.408 nm)	≤ 36.00	25.74	
Ba (614.171 nm)	≤ 42.00	25.23	
Ar (675.283 nm)	≤ 74.00	58.92	
K (766.491 nm)	≤ 80.00	63.16	

Page 2 of 4

เอกสารไม่ควบคุม

Sensitivity Test						Pass
Radial						
Element Wavelength	Specification	Method	Ratio	Standard	Blank	
As (188.980 nm)	≥ 40.0	SRBR	110.5	868.9	54.3	
Se (196.026 nm)	≥ 41.0	SRBR	88.3	934.7	91.3	
Zn (213.857 nm)	≥ 1421.0	SRBR	3535.4	44017.7	153.9	
Pb (220.353 nm)	≥ 46.0	SRBR	184.5	2492.3	159.8	
Mn (257.610 nm)	≥ 3518.0	SRBR	11099.6	249595.3	503.6	
Al (396.152 nm)	≥ 3.4	SBR	8.7	50274.4	5172.0	
Ba (493.408 nm)	≥ 34.0	SBR	124.5	1903164.1	15198.0	
K (766.491 nm)	≥ 1.8	SBR	8.9	110041.4	13891.2	
Axial						
Element Wavelength	Specification	Method	Ratio	Standard	Blank	
As (188.980 nm)	≥ 208.0	SRBR	253.3	3744.3	196.3	
Se (196.026 nm)	≥ 159.0	SRBR	206.7	4199.7	347.2	
Zn (206.200 nm)	≥ 234.0	SRBR	923.0	12282.3	172.1	
Zn (213.857 nm)	≥ 1743.0	SRBR	6368.3	157551.5	601.7	
Cd (214.439 nm)	≥ 4227.0	SRBR	5089.2	99873.7	385.2	
Pb (220.353 nm)	≥ 320.0	SRBR	389.0	10641.1	658.6	
Mn (257.610 nm)	≥ 10625.0	SRBR	21190.4	985528.7	2153.6	
Cr (267.716 nm)	≥ 1048.0	SRBR	3054.1	131797.6	1811.5	
Cu (324.754 nm)	≥ 19.0	SBR	36.3	301401.4	8082.9	
Al (396.152 nm)	≥ 6.0	SBR	10.8	228359.5	19280.5	
Ba (493.408 nm)	≥ 80.0	SBR	106.5	8460421.5	60122.8	
K (766.491 nm)	≥ 24.0	SBR	30.2	1639840.6	52562.1	

Page 3 of 4

เอกสารไม่ควบคุม

Precision Test			Pass
Radial			
Element Wavelength	Specification	Measured Value % RSD	
As (188.980 nm)	≤ 2.60	1.56	
Se (196.026 nm)	≤ 2.60	1.16	
Zn (213.857 nm)	≤ 1.50	0.50	
Pb (220.353 nm)	≤ 2.80	0.74	
Mn (257.610 nm)	≤ 1.50	0.63	
Al (396.152 nm)	≤ 1.50	0.54	
Ba (493.408 nm)	≤ 1.50	0.78	
K (766.491 nm)	≤ 1.50	0.44	
Axial			
Element Wavelength	Specification	Measured Value % RSD	
As (188.980 nm)	≤ 1.50	0.82	
Se (196.026 nm)	≤ 1.50	0.82	
Zn (206.200 nm)	≤ 1.50	0.35	
Zn (213.857 nm)	≤ 1.50	0.34	
Cd (214.439 nm)	≤ 1.50	0.44	
Pb (220.353 nm)	≤ 1.50	0.48	
Mn (257.610 nm)	≤ 1.50	0.83	
Cr (267.716 nm)	≤ 1.50	0.53	
Cu (324.754 nm)	≤ 1.50	0.69	
Al (396.152 nm)	≤ 1.50	0.56	
Ba (493.408 nm)	≤ 1.50	1.29	
K (766.491 nm)	≤ 1.50	0.74	

Page 4 of 4

เอกสารไม่ควบคุม