

ภาคผนวก ง

ใบรับรองการสอบเทียบเครื่องมือ



right solutions.
right partner.

รายการเครื่องมือที่ใช้ในการวิเคราะห์ / ทดสอบ

Sample Name	Parameter	Equipment Name	ID No.	Calibrated Date	Next Cal	Freq. Calibrate (Months)
Ambient	Particulate Matter (PM-10)	High Volume	SGK_FS0122	-	-	On site Calibration
Ambient	Particulate Matter (PM-10)	Digital Balance	BKK_EN0403	3-Jun-24	3-Jun-25	12
Ambient	Total Suspended Particulate	High Volume	SGK_FS0085	-	-	On site Calibration
Ambient	Total Suspended Particulate	Digital Balance	BKK_EN0403	3-Jun-24	3-Jun-25	12
Songkhla Lab	Arsenic	ICP-MS	SGK_CL0048	4-Feb-25	3-Aug-26	18
Songkhla Lab	Arsenic	Cold Room	SGK_CL0065	1-Jul-24	1-Jan-26	18
Songkhla Lab	Cadmium	ICP-MS	SGK_CL0048	4-Feb-25	3-Aug-26	18
Songkhla Lab	Cadmium	Cold Room	SGK_CL0065	1-Jul-24	1-Jan-26	18
Songkhla Lab	Chromium	ICP-MS	SGK_CL0048	4-Feb-25	3-Aug-26	18
Songkhla Lab	Chromium	Cold Room	SGK_CL0065	1-Jul-24	1-Jan-26	18
Songkhla Lab	Copper	ICP-MS	SGK_CL0048	4-Feb-25	3-Aug-26	18
Songkhla Lab	Copper	Cold Room	SGK_CL0065	1-Jul-24	1-Jan-26	18
Songkhla Lab	Iron	ICP-MS	SGK_CL0048	4-Feb-25	3-Aug-26	18
Songkhla Lab	Iron	Cold Room	SGK_CL0065	1-Jul-24	1-Jan-26	18
Songkhla Lab	Lead	ICP-MS	SGK_CL0048	4-Feb-25	3-Aug-26	18
Songkhla Lab	Lead	Cold Room	SGK_CL0065	1-Jul-24	1-Jan-26	18
Songkhla Lab	Manganese	ICP-MS	SGK_CL0048	4-Feb-25	3-Aug-26	18
Songkhla Lab	Manganese	Cold Room	SGK_CL0065	1-Jul-24	1-Jan-26	18
Songkhla Lab	Zinc	ICP-MS	SGK_CL0048	4-Feb-25	3-Aug-26	18
Songkhla Lab	Zinc	Cold Room	SGK_CL0065	1-Jul-24	1-Jan-26	18
Songkhla Lab	Mercury	ICP-MS	SGK_CL0048	4-Feb-25	3-Aug-26	18
Songkhla Lab	Mercury	Cold Room	SGK_CL0065	1-Jul-24	1-Jan-26	18
Songkhla Lab	Total Coliform	Autoclave	SGK_ML0001	1-Jul-24	1-Jan-26	18
Songkhla Lab	Total Coliform	Incubator	SGK_ML0013	15-Jan-24	15-Jul-25	18
Songkhla Lab	Total Coliform	pH Meter	SGK_ML0016	2-Jul-24	2-Jan-26	18
Songkhla Lab	Total Coliform	Water Bath	SGK_ML0021	1-Jul-24	1-Jan-26	18
Songkhla Lab	<i>Escherichia coli</i>	Autoclave	SGK_ML0002	1-Jul-24	1-Jan-26	18
Songkhla Lab	<i>Escherichia coli</i>	Incubator	SGK_ML0013	15-Jan-24	15-Jul-25	18
Songkhla Lab	<i>Escherichia coli</i>	pH meter	SGK_ML0016	2-Jul-24	2-Jan-26	18
Songkhla Lab	<i>Escherichia coli</i>	Water Bath	SGK_ML0021	1-Jul-24	1-Jan-26	18
Water Lab	Fluoride	Ion Chromatography	BKK_EN0427	21-Nov-24	21-Nov-25	12
Water Lab	Nitrate	Ion Chromatography	BKK_EN0427	21-Nov-24	21-Nov-25	12
Water Lab	Chloride	Ion Chromatography	BKK_EN0427	21-Nov-24	21-Nov-25	12
Water Lab	Sulfate	Ion Chromatography	BKK_EN0427	21-Nov-24	21-Nov-25	12
Songkhla Lab	BOD	Incubator	SGK_CL0028	10-Jan-25	10-Jul-26	18
Songkhla Lab	BOD	DO/BOD Analyser	SGK_CL0073	21-May-24	21-Nov-25	18
Songkhla Lab	pH at 25 °C	pH meter	SGK_CL0030	19-Oct-24	19-Apr-26	18
Songkhla Lab	Oil & Grease	Electronic Top-Loading Balance	SGK_CL0045	10-Jan-25	10-Jan-26	12
Songkhla Lab	Oil & Grease	Oven	SGK_CL0024	19-Oct-24	19-Apr-26	18
Songkhla Lab	Oil & Grease	Water Bath	SGK_CL0035	10-Jan-25	10-Jul-26	18
Songkhla Lab	Sulfide	Cold Room	SGK_CL0065	1-Jul-24	1-Jan-26	18
Songkhla Lab	Total Dissolved Solids 180°C	Electronic Top-Loading Balance	SGK_CL0045	10-Jan-25	10-Jan-26	12
Songkhla Lab	Total Dissolved Solids 180°C	Oven	SGK_CL0024	19-Oct-24	19-Apr-26	18
Songkhla Lab	Total Suspended Solids	Electronic Top-Loading Balance	SGK_CL0045	10-Jan-25	10-Jan-26	12
Songkhla Lab	Total Suspended Solids	Oven	SGK_CL0024	19-Oct-24	19-Apr-26	18
Songkhla Lab	Total Kjeldahl Nitrogen	KJELDATHERM@BLOCK-DIGESTION	SGK_CL0092	25-Feb-25	25-Feb-26	12
Songkhla Lab	Fecal Coliform	Autoclave	SGK_ML0001	1-Jul-24	1-Jan-26	18
Songkhla Lab	Fecal Coliform	Incubator	SGK_ML0013	15-Jan-24	15-Jul-25	18
Songkhla Lab	Fecal Coliform	pH Meter	SGK_ML0016	2-Jul-24	2-Jan-26	18
Songkhla Lab	Fecal Coliform	Water Bath	SGK_ML0021	1-Jul-24	1-Jan-26	18
Songkhla Lab	Total Alkalinity	pH meter	SGK_CL0030	19-Oct-24	19-Apr-26	18
Songkhla Lab	<i>Staphylococcus aureus</i>	Autoclave	SGK_ML0002	1-Jul-24	1-Jan-26	18
Songkhla Lab	<i>Staphylococcus aureus</i>	Incubator	SGK_ML0013	15-Jan-24	15-Jul-25	18



right solutions.
right partner.

รายการเครื่องมือที่ใช้ในการวิเคราะห์ / ทดสอบ


Sample Name	Parameter	Equipment Name	ID No.	Calibrated Date	Next Cal	Freq. Calibrate (Months)
Songkhla Lab	<i>Staphylococcus aureus</i>	pH meter	SGK_ML0016	2-Jul-24	2-Jan-26	18
Water Lab	Ammonia Nitrogen	Discrete analyzer	BKK_EN0037	16-Aug-24	16-Aug-25	12
Water Lab	<i>Pseudomonas aeruginosa</i>	Autoclave	BKK_ML0043	26-Sep-24	26-Mar-26	18
Water Lab	<i>Pseudomonas aeruginosa</i>	Incubator	BKK_ML0010	3-Dec-24	3-Dec-25	12
Water Lab	<i>Pseudomonas aeruginosa</i>	Hot Air Oven	BKK_ML0013	23-Apr-24	23-Oct-25	18
Songkhla Lab	Turbidity	Turbidity Meter	SGK_FS0045	22-May-24	22-May-25	12
Water Lab	Cyanuric acid	Spectrophotometer	SGK_CL0100	25-Dec-24	25-Dec-25	12



High Volume Air Sampler Calibration Worksheet

Project Site : Southern Star Hotel Co.,Ltd. Barometric Pressure (mm Hg) : 755.3
Calibrate Location : บริเวณพื้นที่โครงการ Temperature (°C) : 32.4
Calibrate Date : 31-Mar-25 High Volume ID : SGK_FS0122
CalibrationSheet No.: C-310325-SGK_FS0122 High Volume Model : TE-5009X
Calibrator ID: SGK_FS0034 High Volume S/N : 1089
Calibrator Model : TE-5028A Calibrator Slope : 0.98245
Calibrator S/N : 3449 Calibrator Intercept : -0.02325

Test No.	Delta H ₂ O (inch)	Qa (m ³ /min)	I : Chart (CFM)	Linear Regression
1	1.8	0.892	32	Slope : 36.7810 Intercept : -0.5737 Correlation Coefficient : 0.9991
2	2.2	0.983	36	
3	2.8	1.106	40	
4	3.4	1.217	44	
5	4.0	1.318	48	

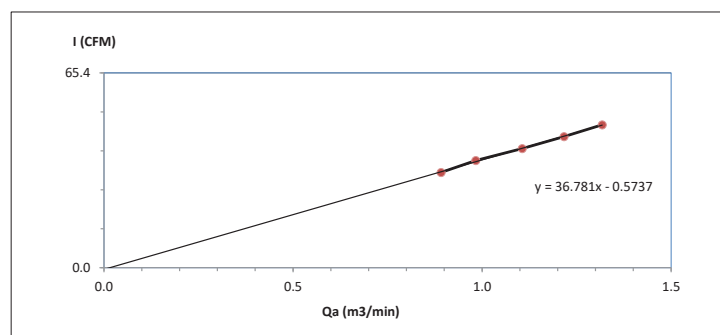

Calibrated by Yong R.
(Mr. Yongsil Rangsee)
Field Services Scientist(2)

Approved by Spt S
(Mr. Supot Salamteh)
Field Services Section Head

FORM NO.: F 06-074 REVISION NO.:2 ISSUE DATE: 20/11/23



PLAY SOLUTION TECHNOLOGY COMPANY LIMITED
179/75 Nawong Pracha Pattana Road, Sikan, Donmuang, Bangkok 10210
Tel.:+66 2 011 0505, Fax:+66 2 010 7700
www.playsotec.com



CERTIFICATE OF CALIBRATION

Certificate No. : PST-0126-24

W/O No. : WO-0051-24

Customer

Company : ALS LABORATORY GROUP (THAILAND) CO., LTD.
Address : 104 Phatthanakan 40, Phatthanakan Road, Khwaeng Phatthanakan,
City / Province : Khet Suan Luang, Bangkok
Zip/Postal : 10250

Page no. 1 of 3

Device

Equipment : Electronic Balance Capacity : 120 / 220 g
Manufacturer : OHAUS Readability : 0.00001 / 0.0001 g
Model : EX225D/AD ID No. : BKK_EN0403
Serial No. : C309774648
Condition : Normal

Environment Conditions

Location of Calibration : Enviroment Lab
Ambient Temperature : 20.1 (°C) ± 3 °C
Relative Humidity : 70.3 (%RH) ± 15 %RH
Barometric Pressure : 1011.1 (mba) ± 10 hPa
Comment :

REVIEW BY Jinda K.
APPROVED BY Siriluk P.
NEXT CAL. DATE 03/06/25

Date of Receipt : June 3, 2024
Date of Calibration : June 3, 2024
Issue Date : June 5, 2024

Calibrated by : Mr.Kittichai Rattanatham
Calibrator

Approved by K. O.
(Mr.Kittichai Rattanatham)
Approved Signature

The reported measurement result relates only to the measurand and applies only at the time of measurement.

This Certificate is issued in accordance with the conditions of accreditation granted by Thai Laboratory Accreditation scheme which has assessed the measurement capability of the laboratory and is traceability to recognize national standards and to the unit of measurement realized at the corresponding national standard laboratory. This certificate may not be reproduced other than in full except with the prior written approval prior written approval of the calibration center, Play Solution Technology Co.,Ltd

F-039

REV.03 30/08/66

CERTIFICATE OF CALIBRATION

Certificate No. : PST-0126-24

W/O No. : WO-0051-24

Result of Calibration : Without Adjustment

Page no. 2 of 3

1. Repeatability

Weighing Range 1	g	Nominal Value	g	Standard Deviation	g
Max.capacity	220	50		0.000012	
		200		0.000048	

2. Linearity, Departure of Indication from nominal value

Weighing Range 1

Nominal Value	Standard Value	Indication	Error of Indication	Expanded Uncertainty	Factor k
g	g	g	g	g	
0.01	0.01000	0.01000	-0.000001	0.000082	2.87
0.1	0.10001	0.10001	0.000004	0.000082	2.87
0.5	0.50000	0.50001	0.000012	0.00008	2.87
1	1.00001	1.00002	0.000013	0.00008	2.87
5	5.00002	5.00003	0.000009	0.00008	2.52
10	9.99999	9.99999	-0.000001	0.00008	2.28
50	50.00001	49.99998	-0.000027	0.00016	2.00
100	100.00002	100.00002	0.000004	0.00030	2.00
150	150.00002	150.00001	0.000077	0.00045	2.00
200	200.00003	200.00001	0.000068	0.00060	2.00

CERTIFICATE OF CALIBRATION

Certificate No. : PST-0126-24

W/O No. : WO-0051-24

Result of Calibration

Page no. 3 of 3

3. Eccentricity

Test load at least 1/3 of the maximum capacity, typically placed between 1/2 and 1/3 of the distance from the centre of the load receptor to the edge.



Weighing Range 1

Test Load : 100 g

Position	Indication	g
1	100.00004	
2	100.00005	
3	100.00002	
4	100.00004	
5	100.00003	
Max.Deviation	0.00002	

Standard method

The calibration was performed by using calibration laboratory's in-house calibration method : CP-M-001 based on "UKAS LAB 14 : Calibration of weighing machine" : edition 6 | October 2019

Reference standards instrument

Instrument	OIML Class	S/N	Certificate No.	Due Date
Standard Weight Set	E2	4000021952	22-128725	November 30, 2024
Standard Weight Set	-	-	-	-
Standard Weight Set	-	-	-	-
Standard Weight Set	-	-	-	-

Measurement Uncertainty

The given measurement uncertainty is the standard of the measurement multiplied by an extension factor k which corresponds to a confidence level of about 95% for a normal distribution. The standard uncertainty was calculated according to UKAS M3003.

Traceability : The measurement is traceable to national standard, which realize the physical unit of measurement (SI)
- Through the reference calibration laboratory of Asia Medical and Agricultural Laboratory and Research Center Co.,Ltd

END OF REPORT



High Volume Air Sampler Calibration Worksheet

Project Site : Southern Star Hotel Co.,Ltd. Barometric Pressure (mm Hg) : 755.3

Calibrate Location : บริเวณพื้นที่โครงการ Temperature (°C) : 32.4

Calibrate Date : 31-Mar-25 High Volume ID : SGK_FS0085


Calibration Sheet No.: C-310325-SGK_FS0085 High Volume Model : TE-5170D

Calibrator ID: SGK_FS0034 High Volume S/N : 5990

Calibrator Model : TE-5028A Calibrator Slope : 1.56856

Calibrator S/N : 3449 Calibrator Intercept : -0.03727

Test No.	Delta H ₂ O (inch)	Q _{std} (m ³ /min)	I : Chart (CFM)	Linear Regression
1	2.8	1.0878	40	Slope : 37.5540 Intercept : -1.0835 Correlation Coefficient : 0.9994
2	3.5	1.2118	44	
3	4.4	1.3542	50	
4	5.2	1.4689	54	
5	6.4	1.6255	60	

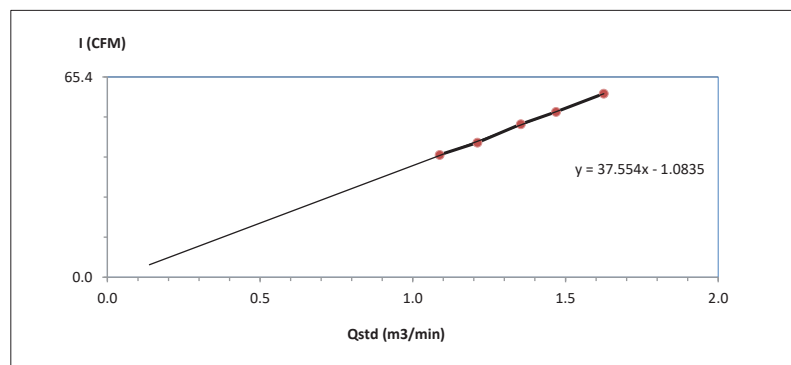

Calibrated by YMR
(Mr. Yongsil Rangsee)
Field Services Scientist(2)

Approved by : Supot S
(Mr. Supot Salamteh)
Field Services Section Head



Agilent Technologies

Agilent Technologies (Thailand) Limited
U CHU LIANG BLDG. 22/F UNIT A,D
968 RAMA 4 ROAD, SILOM, BANGRAK
Bangkok 10500 Thailand

Tel. +662 637 6363
Fax: +662 632 4334
Email: ccc-smt@agilent.com
Website: www.agilent.com/chem

Customer Contact:

ALS Laboratory Group (Thailand) Co
Ltd Head Office

104 Phatthanakan 40 Phatthanakan Rd
Khwaeng Phatthanakan Khet Suan

TAX ID : 0105540004859

tanyatarn.mongkonjirawut@alsglobal.com
027603077

Invoice To:

ALS Laboratory Group (Thailand) Co
Ltd Head Office

104 Phatthanakan 40 Phatthanakan Rd
Khwaeng Phatthanakan Khet Suan

Payer:

ALS Laboratory Group (Thailand) Co
Ltd Branch Number 0002
114/1 Moo8 Banplu Subdistric Hat
Yai District

SONGKHLA 90250

Delivery Site:


ALS Laboratory Group (Thailand) Co
Ltd
Branch Number 0002
114/1 Moo8 Banplu Subdistric Hat
Yai District

Location:

Room
Bldg
Lab
Dept

SERVICE REPORT

Customer Purchase Order Number:	Customer Number: 70371013
Service Request:	Service Request Date:
Service Order: 6007493912	Service Confirmation: 6906126918

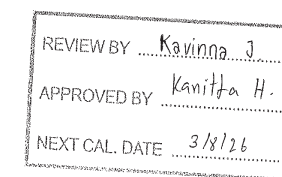

Direct Inquiries to:

Contact Name: Customer Contact Center
Contact E-mail: ccc-smt@agilent.com
Contact Telephone: +662 637 6363
Contact Fax: +662 632 4334

products | applications | software | services

Learn more about Agilent's Special Offers, Products, Services and our
full range of laboratory productivity solutions optimized for your
applications and workflows. Visit us at www.agilent.com/chem

Agilent Technologies (Thailand) Limited. Head Office
U Chu Liang Bldg. 22/F Unit A,D
968 Rama 4 Road, Silom, Bangrak,
Bangkok 10500 Thailand
Tax ID : 0105542068218

Citibank N.A. Bangkok Branch ORIGINAL
399 Interchange 21 Building, Sukhumvit Road, Klongtoey Nau
Sub-district, Wattana District, Bangkok 10110 Thailand
Acc. No: 012-4452-007.
THB:Krung Thai Bank PCL
Siam Square Br.,416/1-2 Rama I Rd.,Pathumwan, BKK 10330
Thailand

Service Confirmation Number: 6906126918
Service Confirmation Date: 04.02.2025

Service Instrument:

Model Number	Model Description	Serial Number	System Handle	Parent Asset
SYS-IM-7900	ICPMS 7900 System			
G7201C	ICP-MS MassHunter SW only (excludes PC)	USH3799575	ICP MS 7900	SYS-IM-7900
G8403A	Agilent 7900 ICP-MS	JP16511669	ICP MS 7900	SYS-IM-7900
G8411A	ISIS 3 for Agilent 7850/7900/8900	JP16510379	ICP MS 7900	SYS-IM-7900


Service Items:

Item	Service/Part #	Description	Qty	Entitlement	Service Start	Service End
1000	EOQ	Enterprise Operational Qualification	1.00	Agreement Entitlement - 100 % covered	03.02.2025	04.02.2025
1010	5185-5850	ICP-MS Checkout Solutions	1.00	Agreement Entitlement - 100 % covered		

Additional Information:

Service Confirmation Number: 6906126918
Service Confirmation Date: 04.02.2025

Service Information:

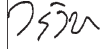Problem Description: *PR-OQ-ICP MS 7900-5001391286		
Service Provided: OQ_SGK_CL0048 Perform oq test . Sensitivity, resolution background both He and H2 mode. Stability test all pass. Sent oq report to customer. Let s customer test run sample.		
Service Overview Code: Reason Code: Scheduled Service Diagnosis Code: Scheduled Service Resolution Code: Scheduled Service		
Reported Hours: 7.0	Travel Hours: 2.0	
Customer Field Service Representative Name: Worawit Timakul	Customer Field Service Representative Signature: 	Date: 04 Feb 2025
Customer Name: TANYATORN MONGKONJIRAWUT	Customer Signature: 	Date: 04 Feb 2025
Additional Comments:		



CARIBRATION CERTIFICATE

Issued Date : 4-Jul-2024

Certificate No. : 24TH2757

CSR No. : A150/07473

Page. : 1 of 3

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Karnchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Chemical Laboratory
Instrument Name : Cold Room
Manufacturer : Danfoss
Model : N/A
Serial No. : N/A
ID No. : SGK_CL0065
Resolution : 0.1 °C
Received Date : 1-Jul-2024
Calibrated Date : 1-Jul-2024
Ambient Temperature : (30 ± 10) °C
Relative Humidity : (50 ± 30) %

REVIEW BY Ananta B..
APPROVED BY Kanitta H..
NEXT CAL DATE 1/01/26

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WI.012 based on GLA - 20

The Southern Calibration Service Co.,Ltd.calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and /or national standards which realize the units of measurement
according to the International System of Unit (SI) through :

- TISTR : Thailand Institute of Scientific and Technological Research

Calibrated by : Ibrorhim Saleemin

Approved by :

Imron Rattanaylum / Technical Manager



The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 24TH2757

CSR No. : A150/07473

Page. : 2 of 3

Details of Calibration


1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Data Acquisition/Switch Unit	34970A	MY58009813	PSL-T0707-1/67	22-May-2025

2. The results reported in this certificate refer to the condition of the instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.
3. This certificate is not certified any commercial transaction
4. Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration : (✓) Without Adjustment () After Adjustment

1. Sensor Installation Diagram

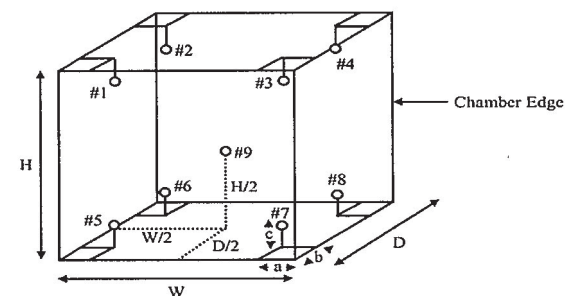

Sensor Installation Details

a = 5.0 cm
b = 5.0 cm
c = 5.0 cm

Dimension of the chamber

W = 40.0 cm
H = 40.0 cm
D = 33.0 cm



Certificate No. : 24TH2757
CSR No. : A150/07473
Page. : 3 of 3

Result of Calibration :

2. Temperature Measurement Accuracy Test

The measurement results of the Cold Room and associates are reported in the manner as shown below

Cal point (°C)	Measured Standard Temperature At Spread Locations (°C)									Uncertainty (± °C)
	#1	#2	#3	#4	#5	#6	#7	#8	Ref. 9	
4	3.90	4.04	3.72	3.78	3.69	3.76	3.55	3.43	3.34	0.39

3. Performance Result

The performance of the Cold Room are reported as shown below

Cal point (°C)	UUC Setting (°C)	UUC Reading (°C)	Temperature Stability (± °C)	Temperature Uniformity (°C)	Overall Variation (°C)
4	4.0	4.0	0.10	0.75	0.75

- UUC = Unit Under Calibration

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$,
providing a level of confidence of approximately 95%

... End ...



Southern Calibration Service Co., Ltd.

669/35 Karnjanavanit Rd., Banpru, Hatyai, Songkla 90250 Thailand

Tel : 08 1599 0417 Fax : 0 7480 5133 Email : s.calibration@gmail.com www.scal-lab.com



NSC-TISI-TIS17025
CALIBRATION 0294

CARIBRATION CERTIFICATE

Issued Date : 4-Jul-2024

Certificate No. : 24TH2760

CSR No. : A150/07473

Page. : 1 of 3

Customer : ALS Laboratory Group (Thailand) Co., Ltd

114/1 Moo 8, Kamchanawanich Rd. Tambon, Ban Phru,

Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Microbiological Laboratory

Instrument Name : Autoclave

Manufacturer : TOMY

Model : SX-700

Serial No. : 52134079

ID No. : SGK_ML0001

Resolution : 1 °C

Received Date : 1-Jul-2024

Calibrated Date : 1-Jul-2024

Ambient Temperature : (30 ± 10) °C

Relative Humidity : (50 ± 30) %

REVIEW BY Wichapol P.

APPROVED BY Kamida H.

NEXT CAL. DATE 1/01/26

Calibration Method Used :
This instrument was calibrated using the Calibration In - house method : SCAL.WI.013 based on BS 2646 : 1993 (part 5)
The Southern Calibration Service Co.,Ltd.calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :
This Certificate is traceable to the International and /or national standards which realize the units of measurement
according to the International System of Unit (SI) through :
- SCaL : Sounthern Calibration Service Co., Ltd.,

Calibrated by : Ibrorhim Saleemin

Approved by : 
Imron Rattanaylum / Technical Manager



The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 24TH2760
CSR No. : A150/07473
Page. : 2 of 3

Details of Calibration


1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Data logger With Sensor	GL240	C90432223	24SDAT005	2-May-2025

2. The results reported in this certificate refer to the condition of the instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.
3. This certificate is not certified any commercial transaction
4. Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration : (✓) Without Adjustment () After Adjustment

1. Sensor Installation Diagram

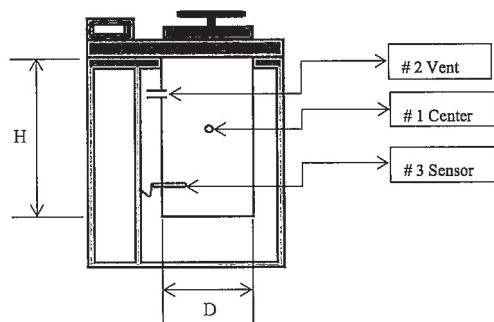

Chamber Diameter (D) : 30 cm

Chamber Height (H) : 70 cm



Certificate No. : 24TH2760
CSR No. : A150/07473
Page. : 3 of 3

Result of Calibration :

2. Temperature Measurement Accuracy Test

The measurement results of the Autoclave and associates are reported in the manner as shown below

Cal point (°C)	Measured Standard Temperature At Spread Locations (°C)			Pressure Reading	Uncertainty (± °C)
	Center #1	Vent #2	Sensor #3		
115	116.3	116.4	116.3	0.07 MPa	0.76
118	119.2	119.2	119.2	0.09 MPa	0.76
121	121.8	121.5	121.6	0.11 MPa	0.76

3. Performance Result

The performance of the Autoclave are reported as shown below

Cal point (°C)	UUC Setting (°C)	UUC Reading (°C)	Temperature Stability (± °C)	Temperature Uniformity (°C)	Overall Variation (°C)
115	115	115	0.10	0.17	0.17
118	118	118	0.10	0.10	0.10
121	121	121	0.70	0.50	0.84

- UUC = Unit Under Calibration

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%

... End ...



CARIBRATION CERTIFICATE

Issued Date : 18-Jan-2024

Certificate No. : 24TH0186

CSR No. : A123/06123

Page. : 1 of 3

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Karnchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Microbiological Laboratory
Instrument Name : Incubator
Manufacturer : Memmert
Model : ICP750
Serial No. : F816.0061
ID No. : SGK_ML0013
Resolution : 0.1 °C
Received Date : 15-Jan-2024
Calibrated Date : 15-Jan-2024
Ambient Temperature : (30 ± 10) °C
Relative Humidity : (50 ± 30) %

REVIEW BY Wichapol R.
APPROVED BY Kampha H.
NEXT CAL. DATE 15/07/25

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WL.012 based on GLA - 20

The Southern Calibration Service Co.,Ltd.calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and/or national standards which realize the units of measurement
according to the International System of Unit (SI) through :

- SCaL : Sounthern Calibration Service Co., Ltd.,

Calibrated by : Ibrorhim Saleemin

Approved by :

Imron Rattanaylum / Technical Manager



The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 24TH0186

CSR No. : A123/06123

Page. : 2 of 3

Details of Calibration


1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Data Acquisiton/Switch Unit	34970A	MY58009813	23SDAT004	23-May-2024

2. The results reported in this certificate refer to the condition of the instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.
3. This certificate is not certified any commercial transaction
4. Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration : (✓) Without Adjustment () After Adjustment

1. Sensor Installation Diagram

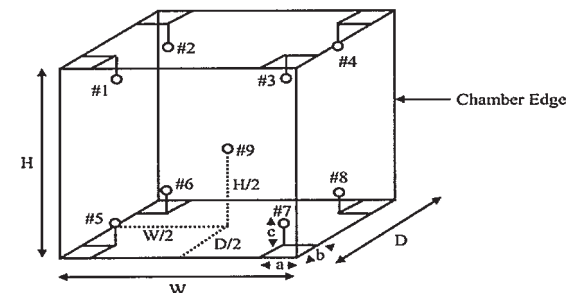

Sensor Installation Details

a = 5.0 cm
b = 5.0 cm
c = 5.0 cm

Dimension of the chamber

W = 104.0 cm
H = 120.0 cm
D = 60.0 cm



Certificate No. : 24TH0186

CSR No. : A123/06123

Page. : 3 of 3

Result of Calibration :**2. Temperature Measurement Accuracy Test**

The measurement results of the Incubator and associates are reported in the manner as shown below

Cal point (°C)	Measured Standard Temperature At Spread Locations (°C)									Uncertainty (± °C)
	#1	#2	#3	#4	#5	#6	#7	#8	Ref. 9	
35	35.04	34.97	34.90	34.89	35.05	34.91	34.99	34.92	34.94	0.38

3. Performance Result

The performance of the Incubator are reported as shown below

Cal point (°C)	UUC Setting (°C)	UUC Reading (°C)	Temperature Stability (± °C)	Temperature Uniformity (°C)	Overall Variation (°C)
35	35.0	35.0	0.10	0.23	0.25

- UUC = Unit Under Calibration

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$,
providing a level of confidence of approximately 95%

... End ...



Southern Calibration Service Co., Ltd.
669/35 Karnjanavanit Rd., Banpru, Hatyai, Songkla 90250 Thailand
Tel : 08 1599 0417 Fax : 0 7480 5133 Email : s.calibration@gmail.com www.scal-lab.com



NAC-MRA
NAC
NSC-TISI-TIS17025
CALIBRATION 0294

CARIBRATION CERTIFICATE

Issued Date : 5-Jul-2024

Certificate No. : 24CH0342

CSR No. : A150/07473

Page. : 1 of 2

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Karnchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Microbiological Laboratory

Instrument Name : pH meter

Manufacturer : Sartorius

Model : PB-10

Serial No. : C07160695

ID No. : SGK_ML0016

Electrode No. : P20067003

Received Date : 2-Jul-2024

Calibrated Date : 2-Jul-2024

Ambient Temperature : (25 ± 3) °C

Relative Humidity : (55 ± 15) %

REVIEW BY Nichaporn P.

APPROVED BY Kanitta H.

NEXT CAL. DATE 2/07/26

Calibration Method Used :
This instrument was calibrated using the Calibration In - house method : SCAL.WI.008 based on direct measurement by
using certified reference Material (CRM)
The Southern Calibration Service Co.,Ltd.calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :
This Certificate is traceable to the International and /or national standards which realize the units of measurement
according to the International System of Unit (SI) through :
- CPAchem : CPAchem Ltd
- WK : WK Electric Co., Ltd.

Calibrated by : Alisara Ma

Approved by : 
Imron Rattanaylum / Technical Manager



The uncertainties are for a confidence probability of approximately 95%
This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 24CH0342
CSR No. : A150/07473
Page. : 2 of 2

Details of Calibration

1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Standard Solution	4.000	61278486	986239	10-May-2025
Standard Solution	7.000	61281486	986241	10-May-2025
Standard Solution	10.01	61287532	986240	10-May-2025
Temperature/Electrical Calibrator	MC2-TE	14987	WK2106-299-223	31-May-2025

2. The results reported in this certificate refer to the condition of the instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.

3. This certificate is not certified any commercial transaction

4. Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration :

1. Electrical Measurement

Applied Voltage (mV)	pH meter Reading		Correction (mV)	Uncertainty (± mV)
	(mV)	(pH)		
177.48	177.5	3.80	-0.02	0.17
0.00	0.0	6.90	0.00	0.13
-177.48	-177.5	10.00	0.02	0.17

2. Before Sample Test Measurement

Standard Buffer Solutions (pH)	pH meter Reading		Correction (pH)	Uncertainty (± pH)
	(pH)	(mV)		
4.008	4.07	164.0	-0.062	0.011
6.985	7.02	-7.1	-0.035	0.020
10.010	10.03	-182.1	-0.020	0.037

3. After Sample Test Measurement

Standard Buffer Solutions (pH)	pH meter Reading		Correction (pH)	Uncertainty (± pH)
	(pH)	(mV)		
4.008	4.02	164.2	-0.012	0.011
6.985	7.00	-0.9	-0.015	0.020
10.010	10.01	-170.7	0.000	0.037

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%

... End ...



Southern Calibration Service Co., Ltd.

669/35 Karnjanavanit Rd., Banpru, Hatyai, Songkla 90250 Thailand

Tel : 08 1599 0417 Fax : 0 7480 5133 Email : s.calibration@gmail.com www.scal-lab.com



CARIBRATION CERTIFICATE

Issued Date : 4-Jul-2024

Certificate No. : 24TH2782

CSR No. : A150/07474

Page. : 1 of 3

Customer : ALS Laboratory Group (Thailand) Co., Ltd

114/1 Moo 8, Karnchanawanich Rd. Tambon, Ban Phru,

Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Microbiological Laboratory

Instrument Name : Liquid Bath

Manufacturer : Memmert

Model : WPE45

Serial No. : L716.0558

ID No. : SGK_ML0021

Resolution : 0.1 °C

Received Date : 1-Jul-2024

Calibrated Date : 1-Jul-2024

Ambient Temperature : (30 ± 10) °C

Relative Humidity : (50 ± 30) %

REVIEW BY Wichapol R.

APPROVED BY Kamitta H.

NEXT CAL. DATE 1/01/26

Calibration Method Used :
This instrument was calibrated using the Calibration In - house method : SCAL.WI.014 based on ASTM E 715 : 1980 (reapproved 2001)
The Southern Calibration Service Co.,Ltd.calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :
This Certificate is traceable to the International and /or national standards which realize the units of measurement according to the International System of Unit (SI) through :
- TISTR : Thailand Institute of Scientific and Technological Research

Calibrated by : Ibrorhim Saleemin

Approved by : 
Imron Rattanaylum / Technical Manager



The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 24TH2782
CSR No. : A150/07474
Page. : 2 of 3

Details of Calibration

1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Data Acquisition/Switch Unit	34970A	MY58009813	PSL-T0707-1/67	22-May-2025


2. The results reported in this certificate refer to the condition of the instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.
3. This certificate is not certified any commercial transaction
4. Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration :

(✓) Without Adjustment

() After Adjustment

1. Sensor Installation Diagram

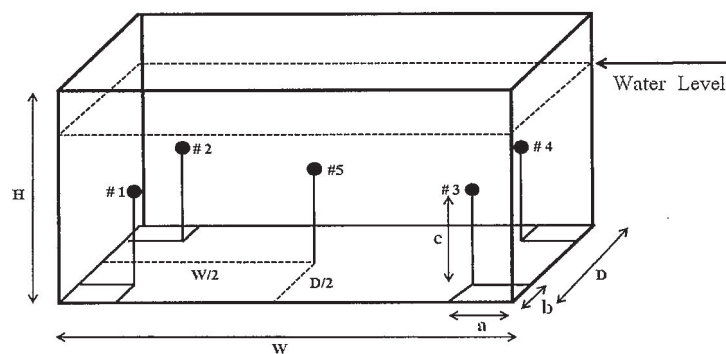

Sensor Installation Details

a = 5 cm
b = 5 cm
c = 5 cm

Dimension of the chamber

W = 45 cm
H = 30 cm
D = 35 cm



Certificate No. : 24TH2782
CSR No. : A150/07474
Page. : 3 of 3

Result of Calibration :

2. Temperature Measurement Accuracy Test

The measurement results of the Liquid Bath and associates are reported in the manner as shown below

Cal point (°C)	Measured Standard Temperature At Spread Locations (°C)					Uncertainty (± °C)
	#1	#2	#3	#4	Ref.5	
44.5	44.48	44.46	44.51	44.49	44.47	0.14

3. Performance Result

The performance of the Liquid Bath are reported as shown below

Cal point (°C)	UUC Setting (°C)	UUC Reading (°C)	Temperature Stability (± °C)	Temperature Uniformity (°C)	Overall Variation (°C)
44.5	44.5	44.5	0.10	0.14	0.25

- UUC = Unit Under Calibration

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%

... End ...



Southern Calibration Service Co., Ltd.

669/35 Karnjanavanit Rd., Banpru, Hatyai, Songkla 90250 Thailand
Tel : 08 1599 0417 Fax : 0 7480 5133 Email : s.calibration@gmail.com www.scal-lab.com



CALIBRATION CERTIFICATE

Issued Date : 4-Jul-2024

Certificate No. : 24TH2761

CSR No. : A150/07473

Page. : 1 of 3

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Karnchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Microbiological Laboratory

Instrument Name : Autoclave

Manufacturer : TOMY

Model : SX-700

Serial No. : 52134080

ID No. : SGK_ML0002

Resolution : 1 °C

Received Date : 1-Jul-2024

Calibrated Date : 1-Jul-2024

Ambient Temperature : (30 ± 10) °C

Relative Humidity : (50 ± 30) %

REVIEW BY Wichapol P.
APPROVED BY Kanitta H.
NEXT CAL. DATE 1/01/26

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WI.013 based on BS 2646 : 1993 (part 5)

The Southern Calibration Service Co.,Ltd.calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and/or national standards which realize the units of measurement according to the International System of Unit (SI) through :

- SCaL : Sounthern Calibration Service Co., Ltd.,

Calibrated by : ibrorhim Saleemin

Approved by :

Imron Rattanaylum / Technical Manager



The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 24TH2761

CSR No. : A150/07473

Page. : 2 of 3

Details of Calibration


1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Data logger With Sensor	GL240	C90432223	24SDAT005	2-May-2025

- The results reported in this certificate refer to the condition of the instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.
- This certificate is not certified any commercial transaction
- Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration : (✓) Without Adjustment () After Adjustment

1. Sensor Installation Diagram

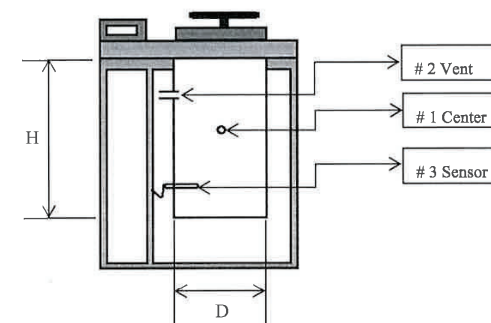

Chamber Diameter (D) : 30 cm

Chamber Height (H) : 70 cm



Certificate No. : 24TH2761
CSR No. : A150/07473
Page. : 3 of 3

Result of Calibration :

2. Temperature Measurement Accuracy Test

The measurement results of the Autoclave and associates are reported in the manner as shown below

Cal point (°C)	Measured Standard Temperature At Spread Locations (°C)			Pressure Reading	Uncertainty (± °C)
	Center #1	Vent #2	Sensor #3		
115	116.3	115.9	116.2	0.07 MPa	0.76
118	119.6	119.4	119.2	0.09 MPa	0.76
121	122.0	122.0	122.1	0.11 MPa	0.76

3. Performance Result

The performance of the Autoclave are reported as shown below

Cal point (°C)	UUC Setting (°C)	UUC Reading (°C)	Temperature Stability (± °C)	Temperature Uniformity (°C)	Overall Variation (°C)
115	115	115	0.10	0.33	0.43
118	118	118	0.30	0.50	0.50
121	121	121	0.20	0.14	0.34

- Operating Time = 900.00 sec

- UUC = Unit Under Calibration

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$,
providing a level of confidence of approximately 95%

... End ...

REVIEW BY	Tanyatorn M.
APPROVED BY	KL AL
NEXT CAL. DATE	21 Nov 25



Bara Scientific
SOLUTION OF SUCCESS

Certificate of Calibration

Shimadzu LC-HIC

This certificate is to verify that instrument below are calibrated
by Bara Scientific Co., Ltd

Instrument	Serial No.
DGU-403	L22166050657
SIL-20AC	L20176012374
CTO-40S	L22236003442
LC-20ADSP	L20106096217
SPD-40	L22256002616
CDD-10Avp	C21346004484
CBM-40lite	L22126103139

For

ALS Laboratory Group (Thailand) Co., Ltd.

Operator Signature: Thawatchai Date: 21 November 2024

(Mr.Thawatchai Toros)

Service Engineer



Bara Scientific Co.,Ltd.
U Chu Liang Building Floor 7 968 Rama 4 Road Silom Bangrak Bangkok 10500
Thailand Tel : 02-6324300 (auto 20 lines) Fax : 02-6375496-7 www.barascientific.com



Southern Calibration Service Co., Ltd.

669/35 Karnjanavanit Rd., Banpru, Hatyai, Songkla 90250 Thailand
Tel : 08 1599 0417 Fax : 0 7480 5133 Email : s.calibration@gmail.com www.scal-lab.com



CALIBRATION CERTIFICATE

Issued Date : 13-Jan-2025

Certificate No. : 25TH0226

CSR No. : A101/05028

Page. : 1 of 3

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Karnchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Customer Laboratory

Instrument Name : Incubator

Manufacturer : Memmert

Model : ICP750

Serial No. : F816.0063

ID No. : SGK_CL0028

Resolution : 0.1 °C

Received Date : 10-Jan-2025

Calibrated Date : 10-Jan-2025

Ambient Temperature : (30 ± 10) °C

Relative Humidity : (50 ± 30) %

REVIEW BY Ananta B.
APPROVED BY Kanida H.
NEXT CAL. DATE 10/1/26

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WI.012 based on GLA - 20

The Southern Calibration Service Co.,Ltd.calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and /or national standards which realize the units of measurement according to the International System of Unit (SI) through :

- TISTR : Thailand Institute of Scientific and Technological Research

Calibrated by : Ibrorhim Saleemin

Approved by :

Imron Rattanaylum / Technical Manager



The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 25TH0226

CSR No. : A101/05028

Page. : 2 of 3

Details of Calibration


1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Data Acquisition/Switch Unit	34970A	MY58009813	PSL-T0707-1/67	22-May-2025

- The results reported in this certificate refer to the condition of the instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.
- This certificate is not certified any commercial transaction
- Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration : (✓) Without Adjustment () After Adjustment

1. Sensor Installation Diagram

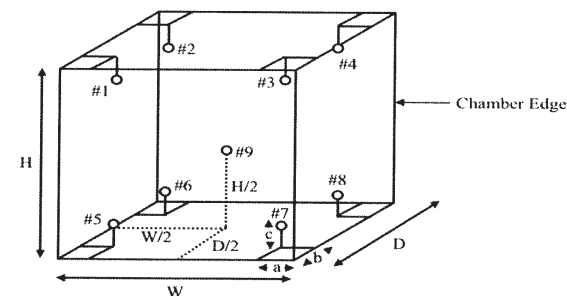

Sensor Installation Details

a = 5.0 cm
b = 5.0 cm
c = 5.0 cm

Dimension of the chamber

W = 104.0 cm
H = 120.0 cm
D = 60.0 cm



Certificate No. : 25TH0226
CSR No. : A101/05028
Page. : 3 of 3

Result of Calibration :

2. Temperature Measurement Accuracy Test

The measurement results of the Incubator and associates are reported in the manner as shown below

Cal point (°C)	Measured Standard Temperature At Spread Locations (°C)									Uncertainty (± °C)
	#1	#2	#3	#4	#5	#6	#7	#8	Ref. 9	
20	20.03	20.23	20.37	20.34	20.32	20.03	20.19	20.31	20.06	0.38

3. Performance Result

The performance of the Incubator are reported as shown below

Cal point (°C)	UUC Setting (°C)	UUC Reading (°C)	Temperature Stability (± °C)	Temperature Uniformity (°C)	Overall Variation (°C)
20	20.0	20.0	0.10	0.31	0.40

- UUC = Unit Under Calibration

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%

... End ...



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3 : EQUIPMENT CALIBRATION AND TESTING SERVICES

534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000 FAX. 0-2719-9484

Certificate of Testing

Cert.No.: 24TW96
Page.: 1 of 2

Equipment : DO Meter
Manufacturer : YSI
Model : 5000
Serial No. : 17B101473
ID No. : SGK_CL0073
Received Date : 17 May 2024
Test Date : 21 May 2024
Reference : 2405-0608DSC-1
Submitted by : ALS Laboratory Group (Thailand) Co.,Ltd.
Songkhla Branch.
114/1 Moo 8 Kamchanawanich Rd., T.Ban Phru,
A.Hat Yai, Songkhla 90250 Thailand
Laboratory Condition : Temperature (25 ± 5) °C
Humidity (50 ± 20) %
Test Procedure : In - house method : CP-CH9
by Comparison Technique with Azide Modification Method
Tested by : Walalak Sirithean
Approved by : 
Approved Signatory
() Unnopphol Harachai
() Ponpan Paipim
(✓) Saithip Meangmai
Issue Date : 21 May 2024

REVIEW BY	Ananta B
APPROVED BY	Kanitta H
NEXT CAL DATE	21/11/25



Cert.No.: 24TW96

Page.: 2 of 2

Condition of this result of calibration

1. Reference Standard Instruments :

This certification is traceable to the International System of Unit through the reference standards laboratory of Industrial Calibration Center, Technology Promotion Association (Thailand-Japan).

Instruments	Serial No.	ID No.	Certificate No.	Due Date
1. Burette	-	130BU10	23CG1172	22 Mar 2025
2. Balance	14233821	110RC001	23MM405	16 July 2024

2. Standard Material :-

Material	Manufacturer	Lot.No.	Assay
Sodium Thiosulfate pentahydrate	Merck	AM1763316	100.2%

Result : Dissolved Oxygen Meter Adjustment With Air 100 %
Dissolved Oxygen Probe No.: 17B100103

Titration Method (Azide Modification Method)	DO Meter Reading	Standard Deviation
(mg/L)	(mg/L)	(mg/L)
8.18	8.18	0.0071

This report was certified only for the instrument we tested. It is allowable to use for study
Intend to use for advertising and referral purpose is prohibited. This report may not be reproduced
other in full, without written approval of the laboratory

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert. No.: 24LM77

Page.: 1 of 2

Equipment : DO Meter with Sensor
Manufacturer : YSI
Model : 5000-115
Serial No. : 17B101473
ID No. : SGK_CL0073
Submitted by : ALS Laboratory Group (Thailand) Co.,Ltd.
Songkhla Branch,
114/1 Moo 8 Kamchanawanich Rd.,
T.Ban Phru, A.Hat Yai,
Songkhla 90250 Thailand
Location : TPA On Site Calibration Laboratory
Received Order : 17 May 2024
Calibrated Date : 27 May 2024
Ambient Temperature : (26 ± 10) °C
Relative Humidity : (50 ± 30) %
AC Line Voltage : (220 ± 22) V

Calibrated by : Khit Ruttanaprapachai

Approved by :

Approved Signatory

() Ponpan Paipim
() Suwit Imjai
(✓) Kunchit Promprat

Issue Date : 28 May 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Equipment : DO Meter with Sensor
Condition As-Received : Used Item
Reference : 2405-0608DSC-2

Cert. No.: 24LM77
Page.: 2 of 2

Procedure Used :-

Calibration were conducted using in-house calibration procedure CP-OT01 according to comparison with Industrial Platinum Resistance Thermometer (IPRT) into Temperature Bath.

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Digital Thermometer	2188080	2311216	TPA	11 Oct 2024

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function : Temperature measurement.

This instrument was connected with temperature sensor, S/N.: 17B100103

Calibration Point	Immersion Depth	Standard Temperature	UUC* Reading	Error	Uncertainty	Coverage Factor
(°C)	(mm)	(°C)	(°C)	(°C)	(± °C)	k
20.00	60	20.005	19.79	-0.215	0.15	2.00

UUC* : Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-



Southern Calibration Service Co., Ltd.

669/35 Karnjanavanit Rd., Banpru, Hatyai, Songkla 90250 Thailand
Tel : 08 1599 0417 Fax : 0 7480 5133 Email : s.calibration@gmail.com www.scal-lab.com



CARIBRATION CERTIFICATE

Issued Date : 22-Oct-2024

Certificate No. : 24CH0526

CSR No. : A163/08133

Page. : 1 of 2

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Karnchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Chemical Laboratory

Instrument Name : pH meter

Manufacturer : Mettler Toledo

Model : S220

Serial No. : B625631849

ID No. : SGK_CL0030

Electrode No. : 2281592

Received Date : 19-Oct-2024

Calibrated Date : 19-Oct-2024

Ambient Temperature : (25 ± 3) °C

Relative Humidity : (55 ± 15) %

REVIEW BY Ananta B.
APPROVED BY Kanitta H.
NEXT CAL. DATE 19/4/26

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WI.008 based on direct measurement by using certified reference Material (CRM)

The Southern Calibration Service Co.,Ltd.calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and/or national standards which realize the units of measurement according to the International System of Unit (SI) through :

- CPAchem : CPAchem Ltd
- WK : WK Electric Co., Ltd.
- SCal : Sounthern Calibration Service Co., Ltd.,

Calibrated by : Alisara Ma

Approved by :

Imron Rattanaylum / Technical Manager



The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 24CH0526

CSR No. : A163/08133

Page. : 2 of 2

Details of Calibration

1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Standard Solution	4.000	61310674	1042701	26-Oct-2025
Standard Solution	7.000	61314184	1042700	26-Oct-2025
Standard Solution	10.01	61313804	1042702	26-Oct-2025
Temperature/Electrical Calibrator	MC2-TE	14987	WK2106-299-223	31-May-2025
Digital Thermometer With Sensor	DP-77	1.360896	24SDTH005	7-Aug-2025

2. The results reported in this certificate refer to the condition of the instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.

3. This certificate is not certified any commercial transaction

4. Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration :

1. Electrical Measurement

Applied Voltage (mV)	pH meter Reading		Correction (mV)	Uncertainty (± mV)
	(mV)	(pH)		
177.48	177.5	4.00	-0.02	0.17
0.00	0.0	7.00	0.00	0.13
-177.48	-177.5	10.00	0.02	0.17

2. Before Sample Test Measurement

Standard Buffer Solutions (pH)	pH meter Reading		Correction (pH)	Uncertainty (± pH)
	(pH)	(mV)		
4.007	3.99	178.1	0.017	0.0092
6.976	7.02	1.7	-0.044	0.019
10.009	9.96	-169.6	0.049	0.038

3. After Sample Test Measurement

Standard Buffer Solutions (pH)	pH meter Reading		Correction (pH)	Uncertainty (± pH)
	(pH)	(mV)		
4.007	3.99	177.7	0.017	0.0092
6.976	7.01	3.7	-0.034	0.019
10.009	10.00	-169.0	0.009	0.038

4. Temperature Measurement

Cal Point (°C)	Standard Temperature (°C)	UUC Reading (°C)	Correction (°C)	Uncertainty (± °C)
25	25.021	25.0	0.02	0.060

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%

... End ...



Southern Calibration Service Co., Ltd.

669/35 Karnjanavanit Rd., Banpru, Hatyai, Songkla 90250 Thailand
Tel : 08 1599 0417 Fax : 0 7480 5133 Email : s.calibration@gmail.com www.scal-lab.com



CALIBRATION CERTIFICATE

Issued Date : 13-Jan-2025

Certificate No. : 25MA0125

CSR No. : A101/05028

Page. : 1 of 3

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Karnchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Customer Laboratory

Instrument Name : Analytical Balance

Manufacturer : Sartorius

Model : MSE224S-100-DU

Serial No. : 34705158

ID No. : SGK_CL0045

Resolution : 0.0001 g

Received Date : 10-Jan-2025

Calibrated Date : 10-Jan-2025

Ambient Temperature : (30 ± 10) °C

Relative Humidity : (50 ± 20) %

REVIEW BY Ananta B
APPROVED BY Kamilla H
NEXT CAL. DATE 10/1/26

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WI.001 based on UKAS LAB 14 : 2015

The Southern Calibration Service Co.,Ltd.calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and /or national standards which realize the units of measurement according to the International System of Unit (SI) through :

- SCaL : Sounthern Calibration Service Co., Ltd.,

Calibrated by : Hadbordee Dettawee

Approved by :

Imron Rattanaylum / Technical Manager



The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 25MA0125
CSR No. : A101/05028
Page. : 2 of 3

Details of Calibration

1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Standard Weight Set	2 mg - 1 kg	11119514/01	24SWS001	3-Jul-2025

2. The results reported in this certificate refer to the condition of the instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.
3. This certificate is not certified any commercial transaction
4. Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration : (✓) Without Adjustment () After Adjustment

1. Repeatability

Nominal Value (g)	Standard Deviation (g)
10	0.00000
20	0.00000
200	0.00000

2. Effect of tare

Nominal Value (g)	Standard Value (g)	Balance Reading (g)	Correction (g)
20	20.0000	20.0000	0.0000
40	40.0001	40.0000	0.0001
60	60.0000	60.0000	0.0000
80	80.0001	80.0000	0.0001
100	100.0000	100.0000	0.0000



Certificate No. : 25MA0125
CSR No. : A101/05028
Page. : 3 of 3


Result of Calibration :

3. Off-centre loading

A mass approximately 100g was placed on a pan and moved to various position .

The balance reading obtained are given in the table.

Position					Maximum Difference (g)
1	2	3	4	5	
100.0000	100.0000	100.0000	100.0000	100.0001	0.0001

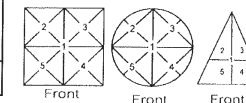

4. Departure from nominal value

Nominal Value (g)	Standard Value (g)	UUC Reading (g)	Correction (g)	Uncertainty (± g)	Coverage Factor (k)
0.01	0.0100	0.0100	0.0000	0.00008	2.0
0.1	0.1000	0.1000	0.0000	0.00008	2.0
0.5	0.5000	0.5000	0.0000	0.00008	2.0
1	1.0000	1.0000	0.0000	0.00008	2.0
2	2.0000	2.0000	0.0000	0.00008	2.0
5	5.0000	5.0000	0.0000	0.00009	2.0
10	10.0000	10.0000	0.0000	0.00009	2.0
20	20.0000	20.0000	0.0000	0.00009	2.0
40	40.0001	40.0000	0.0001	0.00011	2.0
50	50.0000	50.0000	0.0000	0.00013	2.0
60	60.0000	60.0000	0.0000	0.00013	2.0
80	80.0000	80.0000	0.0000	0.00016	2.0
100	100.0000	100.0001	-0.0001	0.00017	2.0
120	120.0000	120.0001	-0.0001	0.00024	2.0
140	140.0001	140.0000	0.0001	0.00024	2.0
160	160.0000	160.0000	0.0000	0.00026	2.0
180	180.0000	180.0000	0.0000	0.00029	2.0
200	200.0000	200.0000	0.0000	0.00030	2.0

- UUC = Unit Under Calibration

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95%

...End...



CARIBRATION CERTIFICATE

Issued Date : 22-Oct-2024

Certificate No. : 24TH4295

CSR No. : A163/08133

Page. : 1 of 3

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Karnchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Chemical Laboratory
Instrument Name : Hot Air Oven
Manufacturer : Memmert
Model : UF110
Serial No. : B416.3387
ID No. : SGK_CL0024
Resolution : 0.1 °C
Received Date : 19-Oct-2024
Calibrated Date : 19-Oct-2024
Ambient Temperature : (30 ± 10) °C
Relative Humidity : (50 ± 30) %

REVIEW BY Anan B.
APPROVED BY Kanitta H.
NEXT CAL. DATE 19/04/2026

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WI.012 based on GLA - 20

The Southern Calibration Service Co.,Ltd.calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and /or national standards which realize the units of measurement
according to the International System of Unit (SI) through :

- TISTR : Thailand Institute of Scientific and Technological Research

Calibrated by : Ibrorhim Saleemin

Approved by :

Imron Rattanaylum / Technical Manager



The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 24TH4295

CSR No. : A163/08133

Page. : 2 of 3

Details of Calibration

1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Data Acquisition/Switch Unit	34970A	MY58009813	PSL-T0707-1/67	22-May-2025


- The results reported in this certificate refer to the condition of the instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.
- This certificate is not certified any commercial transaction
- Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration :

(✓) Without Adjustment

() After Adjustment

1. Sensor Installation Diagram

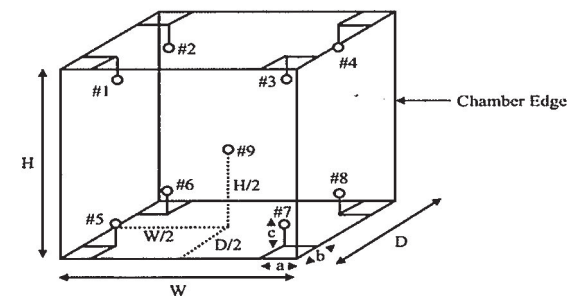

Sensor Installation Details

a = 5.0 cm
b = 5.0 cm
c = 5.0 cm

Dimension of the chamber

W = 55.0 cm
H = 48.0 cm
D = 40.0 cm



Certificate No. : 24TH4295
CSR No. : A163/08133
Page. : 3 of 3

Result of Calibration :

2. Temperature Measurement Accuracy Test

The measurement results of the Hot Air Oven and associates are reported in the manner as shown below

Cal point (°C)	Measured Standard Temperature At Spread Locations (°C)									Uncertainty (± °C)
	#1	#2	#3	#4	#5	#6	#7	#8	Ref. 9	
40	40.36	40.40	40.51	40.43	40.05	40.24	40.09	40.14	39.75	0.38
70	70.27	70.30	70.45	70.24	70.24	70.43	70.29	70.30	69.95	0.36
103	102.94	102.90	103.55	102.96	103.22	103.14	103.10	103.01	102.88	0.36
104	104.15	103.99	104.27	104.06	104.09	104.23	104.26	104.15	103.90	0.36
105	105.04	104.90	105.05	104.87	104.91	104.80	104.82	104.98	104.70	0.36
180	179.19	178.93	179.82	179.10	179.27	179.68	179.12	179.73	179.12	0.41

3. Performance Result

The performance of the Hot Air Oven are reported as shown below

Cal point (°C)	UUC Setting (°C)	UUC Reading (°C)	Temperature Stability (± °C)	Temperature Uniformity (°C)	Overall Variation (°C)
40	40.0	40.0	0.20	0.84	0.84
70	70.0	70.0	0.10	0.59	0.59
103	103.0	103.0	0.20	0.73	0.74
104	104.0	104.0	0.20	0.47	0.56
105	105.0	105.0	0.20	0.44	0.46
180	180.0	180.0	0.50	0.86	1.11

- UUC = Unit Under Calibration

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$,
providing a level of confidence of approximately 95%

... End ...



Southern Calibration Service Co., Ltd.

669/35 Karnjanavanit Rd., Banpru, Hatyai, Songkla 90250 Thailand
Tel : 08 1599 0417 Fax : 0 7480 5133 Email : s.calibration@gmail.com www.scal-lab.com



CALIBRATION CERTIFICATE

Issued Date : 13-Jan-2025

Certificate No. : 25TH0227

CSR No. : A101/05028

Page. : 1 of 3

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Karnchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Customer Laboratory

Instrument Name : Liquid Bath

Manufacturer : Memmert

Model : WNE29

Serial No. : L616.0538

ID No. : SGK_CL0035

Resolution : 0.1 °C

Received Date : 10-Jan-2025

Calibrated Date : 10-Jan-2025

Ambient Temperature : (30 ± 10) °C

Relative Humidity : (50 ± 30) %

REVIEW BYAnanta B.....
APPROVED BYKamilla H.....
NEXT CAL. DATE 10/1/26

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SCAL.WI.014 based on ASTM E 715 : 1980
(reapproved 2001)

The Southern Calibration Service Co.,Ltd.calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and /or national standards which realize the units of measurement
according to the International System of Unit (SI) through :

- TISTR : Thailand Institute of Scientific and Technological Research

Calibrated by : Ibrorhim Saleemin

Approved by :

Imron Rattanaylum / Technical Manager



The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 25TH0227
CSR No. : A101/05028
Page. : 2 of 3

Details of Calibration


1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Data Acquisition/Switch Unit	34970A	MY58009813	PSL-T0707-1/67	22-May-2025

2. The results reported in this certificate refer to the condition of the instrument on the date of calibration and carry no implication regarding the longterm stability of instrument.
3. This certificate is not certified any commercial transaction
4. Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration : (✓) Without Adjustment () After Adjustment

1. Sensor Installation Diagram

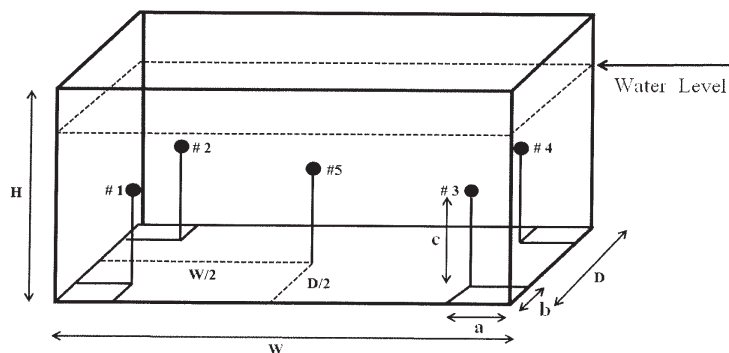

Sensor Installation Details

a = 5 cm
b = 5 cm
c = 5 cm

Dimension of the chamber

W = 59.0 cm
H = 14.0 cm
D = 35.0 cm



Certificate No. : 25TH0227
CSR No. : A101/05028
Page. : 3 of 3

Result of Calibration :

2. Temperature Measurement Accuracy Test

The measurement results of the Liquid Bath and associates are reported in the manner as shown below

Cal point (°C)	Measured Standard Temperature At Spread Locations (°C)					Uncertainty (± °C)
	#1	#2	#3	#4	Ref.5	
80	80.05	80.07	79.98	80.05	80.06	0.17

3. Performance Result

The performance of the Liquid Bath are reported as shown below

Cal point (°C)	UUC Setting (°C)	UUC Reading (°C)	Temperature Stability (± °C)	Temperature Uniformity (°C)	Overall Variation (°C)
80	80.0	80.0	0.20	0.17	0.26

- UUC = Unit Under Calibration

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k = 2, providing a level of confidence of approximately 95%

... End ...



Southern Calibration Service Co., Ltd.

669/35 Karnjanavanit Rd., Banpru, Hatyai, Songkla 90250 Thailand
Tel : 08 1599 0417 Fax : 0 7480 5133 Email : s.calibration@gmail.com www.scal-lab.com

CALIBRATION CERTIFICATE

Issued Date : 28-Feb-2025

Certificate No. : 25TH1073

CSR No. : A185/09219

Page. : 1 of 3

Customer : ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Karnchanawanich Rd. Tambon, Ban Phru,
Amphoe Hat Yai, Songkhla, 90250

Calibration Place : Customer Laboratory

Instrument Name : Digester Block

Manufacturer : Gerhardt

Model : KT-L205BS

Serial No. : 5770230038


ID No. : SGK_CL0092

Receipt Date : 25-Feb-2025

Calibration Date : 25-Feb-2025

Ambient Temperature : (30 ± 10) °C

Relative Humidity : (50 ± 30) %

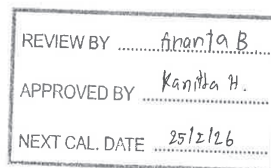

Calibration Method Used :

This instrument was calibrated using the Calibration In - house method : SC.WI.12.016 based on G-20

The Southern Calibration Service Co.,Ltd.calibration control system complies with requirement of ISO/IEC 17025:2017

Traceability of measurement :

This Certificate is traceable to the International and /or national standards which realize the units of measurement according to the International System of Unit (SI) through :

- SCaL : Sounthern Calibration Service Co., Ltd.,

Calibrated by : Ibrorhim Saleemin

Approved by :

Imron Rattanaylum / Technical Manager



The uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of Southern Calibration Service Co., Ltd.



Certificate No. : 25TH1073

CSR No. : A185/09219

Page. : 2 of 3

Details of Calibration

1. Reference Standard Equipment Used:

Equipment	Model	Serial No.	Cert. no.	Due Date
Data Acquisition/Switch Unit	34970A	MY58009813	PSL-T0707-1/67	22-May-2025

2. The results reported in this certificate refer to the condition of the instrument on the date of calibration


and carry no implication regarding the longterm stability of instrument.

3. This certificate is not certified any commercial transaction

4. Condition of Item : normal condition , no indication for any damage or malfunction

Result of Calibration :

Sensor Installation Diagram

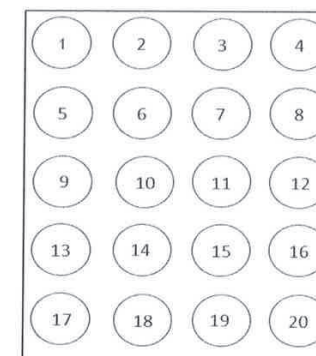

Measurement for Heating Block

Position	UUC Reading (°C)	Standard Temperature (°C)	Correction (°C)	Uncertainty (± °C)
No.1	383	380.24	-2.76	0.32
No.2	383	379.87	-3.13	0.32
No.3	383	379.92	-3.08	0.32
No.4	383	379.94	-3.06	0.32
No.5	383	379.88	-3.12	0.32



Certificate No. : 25TH1073

CSR No. : A185/09219

Page. : 3 of 3

Result of Calibration : (Cont.)

Measurement for Heating Block

Position	UUC Reading (°C)	Standard Temperature (°C)	Correction (°C)	Uncertainty (± °C)
No.6	383	379.97	-3.03	0.32
No.7	383	380.42	-2.58	0.32
No.8	383	380.47	-2.53	0.32
No.9	383	380.21	-2.79	0.32
No.10	383	380.44	-2.56	0.32
No.11	383	379.24	-3.76	0.32
No.12	383	380.27	-2.73	0.32
No.13	383	380.28	-2.72	0.32
No.14	383	379.93	-3.07	0.32
No.15	383	379.92	-3.08	0.32
No.16	383	379.98	-3.02	0.32
No.17	383	379.84	-3.16	0.32
No.18	383	379.82	-3.18	0.32
No.19	383	379.81	-3.19	0.32
No.20	383	380.21	-2.79	0.32

Remark : Calibration Points : 380 °C

UUC Setting : 383 °C

- UUC = Unit Under Calibration

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k = 2,
providing a level of confidence of approximately 95%

...End...

บริษัท ดับเบิล เอส ไดแอกโนสติกส์ จำกัด
DOUBLE S DIAGNOSTICS CO., LTD.4 ซอยอุดมสุข 14 แขวงบางนา เขตบางนา กรุงเทพฯ 10260 โทรศัพท์: (02) 747-7009 โทรสาร: (02) 747-7008
4 Soi Udomsak 14, Bangna, Bangkok 10260 Tel. (02) 747-7009 Fax: (02) 747-7008

Maintenance Plan YEAR : 2024

เดือน	Jan	Feb	Mar	Apr	May	Jun	Jul	Aug	Sep	Oct	Nov	Dec
รวม								5				

Periodical maintenance check list for Konelab

	6M	12M	Note!
1.Diluent-wash tubing change	☐	☒	
2.ISE tubing change	☒	☒	None
3.Syringe check/change		☒	
4.Dispensing check/ change		☒	
5.Waste tubing change when necessary		☒	
6.Lamp check/change	☐	☒	
7.Mixer paddle/paddle change(not Konelab20)		☒	
8.ISE needles check/change		☒	None
9.Pump tubing check/ chance	☐	☒	
10.Broken/worn out part check /change		☒	
11.Peristaltic pump check /cleaning/ lubrication	☐	☒	
12.Heating check		☒	
13.Cooling check		☒	
14.Dispenser mechanic check/adjustment	☐	☒	
15.Cuvette transfer mechanic check/adjustment	☐	☒	
16.Dispenser movement check/adjustment	☐	☒	
17.Sample/reagent register check/adjustment	☐	☒	
18.Dispensing tubing tightness check	☐	☒	
19.Photometer and optics cleaning/check/adjustment	☐	☒	
20.Workstation PC cleaning if necessary	☐	☒	
21.Mechanic cleaning/lubrication	☐	☒	
22.Instrument cleaning if necessary	☐	☒	
23.Complete analyzer testing with waterblank/QC or sample	☐	☒	
24.Test parameters/Adjustment/config. Save to USB key	☐	☒	
25.UPS Test	☐	☒	

Place: ALS LAB Instrument: K20 Aquaken
Date/Time: 16/8/62 Serial no: 22781
Service done by: [Signature] Install date:
Signature of customer: [Signature] Date/Time: 16/08/2024

Laboratory
Analyzer User

8/16/2024 14:53

Performed8/16/2024
LotWB34

=====

ACCEPTANCE CRITERIA

=====

	Result	Limit	Warning
Temperature (°C)	37.8	37.0 +/- 1.0	
Dispensing ratio	16.4	14.8 - 17.2	
CV%	0.29	<1.7	
Photometric noise			
Max SD L340_2 (mA)	0.17	<2.0	
Max SD L340_4 (mA)	0.87	<3.0	
Linearity of photometer			
Slope	1.0141	0.94 - 1.06	
Curvature	0.0053	+/- 0.02	
Max bias from linear fit (mA)	4.3	<15.0	
Max delta %	-1.6	+/- 6.0	
Linearity of sample dispensing			
Proport. volume XDISP2 (?l)	2.06	1.96 - 2.16	
Proport. volume XDISP4 (?l)	4.14	3.85 - 4.40	
XDISP2 CV%	1.21	<2.0	
XDISP4 CV%	0.90	<2.0	
XDISP10 CV%	0.68	<2.0	
Needle 0 ?l volume			
Average (A)	0.005	<0.050	
Standard deviation (A)	0.002	<0.005	
Volume (?l)	0.03	<0.32	

=====

OTHER INFORMATION

=====

Dispensing ratio		Photom. noise: SD (mA)	
Posit	Result (A)	Posit	L340_2L340_4
1	0.1549	1	0.150.80
2	0.1549	2	0.170.79
3	0.1537	3	0.040.65
4	0.1547	4	0.160.31
5	0.1547	5	0.110.58
6	0.1545	6	0.140.87

Laboratory
Analyzer User


8/16/2024 14:53

Linearity of sample dispensing

Test	Absorbance (A)
XDISP2	0.306
XDISP4	0.612
XDISP10	1.471

Linearity of photometer

L340_	Target (A)	Meas (A)	Delta (A)	Delta %
1	0.002	0.006	-0.004	-217.7
2	0.486	0.493	-0.007	-1.5
3	1.451	1.469	-0.018	-1.2
4	1.936	1.963	-0.027	-1.4
5	2.415	2.454	-0.039	-1.6

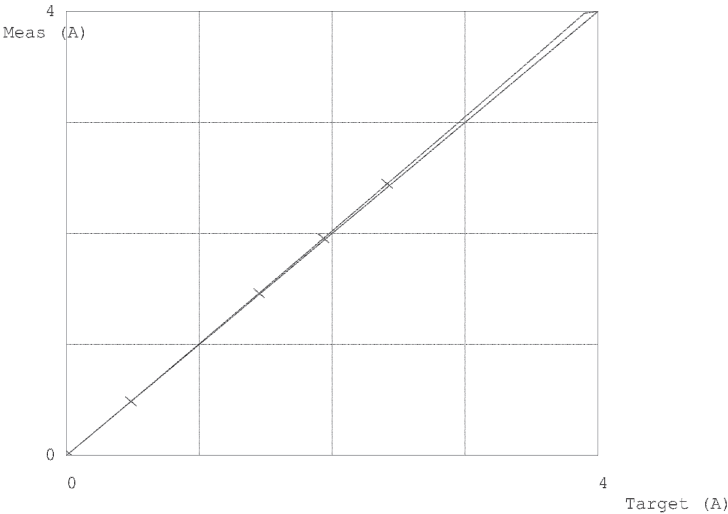



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484

Certificate of Calibration

Cert. No.: 24TM1318
Page : 1 of 3

Equipment : Autoclave
Manufacturer : AES CHEMUNEX
Model : Masterclave 528
Serial No. : 34677152
ID No. : BKK_ML0043

Submitted by : ALS Laboratory Group (Thailand) Co.,Ltd.
104 Phatthanakan 40, Phatthanakan Rd.,
Khwaeng Phatthanakan, Khet Suan Luang,
Bangkok 10250 Thailand

Location : Media Preparation Room

Received Order : 26 September 2024
Calibration Date : 26 September 2024
Ambient Temperature : (26 ± 10) °C
Relative Humidity : (50 ± 30) %

Calibrated by : Tawatchai Pama

Approved by :

Kunchit

Approved Signatory

() Ponpan Paipim
() Suwit Imjai
(✓) Kunchit Promprat

Issue Date : 09 October 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Calibration and Testing Equipment Services.



Equipment : Autoclave
Condition As-Received : Used Item
Reference : 2409-0695OC-6
Procedure Used :-

Cert. No.: 24TM1318
Page : 2 of 3

Calibration were conducted using in-house calibration procedure CP-OT03 Based on BS 2646-5 according to direct measurement method with Data Acquisition which connected with Thermocouple Type T
The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Data Acquisition	MY49001451	24LM79	TPA	29 May 2025

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

4. This result of calibration covers laboratory autoclaves for the sterilization of goods and material which could be infected with organisms categorized as Hazard Group 1, 2 and 3**


(** = Categorization of pathogens according to hazard and categories of containment, second edition, 1990)
It does not cover autoclaves for use with material infect with organisms in Hazard Group 4, for which complete containment and sterilization of infected condensate is considered to be essential.

This result of calibration does not apply to sterilizers or disinfectors used for medical, dental, pharmaceutical or veterinary purposes which are directly concerned with patient care, or those used for fabrics subjected to sterilization which are required to be dry at the end of cycle.

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

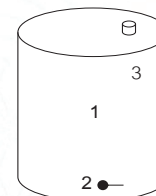

	Environmental		
	(°C)	(%R.H.)	(Volt)
Beginning of Calibration	25	55	225
Finished of Calibration	26	50	226

Position	Description	Ref. Std. ID No.:
1 =	Center of chamber	24-19TC-01
2 =	Temperature sensor	24-19TC-02
3 =	Exhaust port	24-19TC-03



Equipment : Autoclave
Condition As-Received : Used Item
Reference : 2409-0695OC-6
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source

Cert. No.: 24TM1318
Page : 3 of 3

Operating parameter Set : Temperature = 121.0 °C
Sterilization period = 15 minute

UUC* Setting (°C)	UUC* Reading (°C)	Position	Average* Standard Reading (°C)	Stability (± °C)	Pressure Reading (bar)	Uncertainty (± °C)	Coverage Factor <i>k</i>
121.0	120.9	1	121.127	0.19	1.5	0.75	2
		2	121.316				
		3	121.292				

Operating parameter Set : Temperature = 38.0 °C
Sterilization period = 30 minute

UUC* Setting (°C)	UUC* Reading (°C)	Position	Average* Standard Reading (°C)	Stability (± °C)	Pressure Reading (bar)	Uncertainty (± °C)	Coverage Factor <i>k</i>
38.0	38.0	1	37.698	0.17	0.0	0.61	2
		2	37.865				
		3	37.879				

Average* : The average of 30 values in each position.

Stability : One-half of the greatest maximum difference of measured temperature at any one probe.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert. No.: 24TM1398
Page : 1 of 3

Equipment : Incubator
Manufacturer : SHEL-LAB
Model : 1915A
Serial No. : 0200599
ID No. : BKK_ML0010

Submitted by : ALS Laboratory Group (Thailand) Co.,Ltd.
104 Phatthanakan 40, Phatthanakan Rd.,
Khwaeng Phatthanakan, Khet Suan Luang,
Bangkok 10250 Thailand
Location : Incubation & Micrological Reading

Received Order : 03 December 2024
Calibration Date : 03 December 2024
Ambient Temperature : (26 ± 10) °C
Relative Humidity : (50 ± 30) %
AC Line Voltage : (220 ± 22) V

Calibrated by : Kunchit Promprat

Approved by :

Approved Signatory

() Pornthippa Tameyakul
() Ponpan Paipim
(✓) Suwit Imjai

Issue Date : 17 December 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

REVIEW BY Sithichok T.

APPROVED BY

NEXT CAL DATE..... 03/12/25



Equipment : Incubator
Condition As-Received : Used Item
Reference : 2412-0004OC-1

Cert. No.: 24TM1398
Page : 2 of 3

Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD).

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

<u>Instrument</u>	<u>Serial No.</u>	<u>Cert. No.</u>	<u>Traceable</u>	<u>Due Date</u>
1) Data Acquisition	MY49023932	24LM119	TPA	27 Jul 2025

2. This certificate is valid only to the item calibrated on date and place of calibration.


3. This certification is traceable to the International System of Unit.

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Close

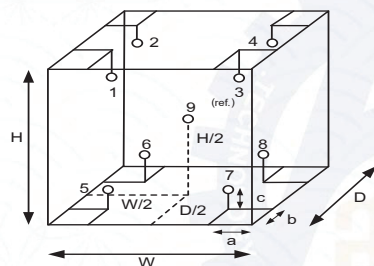

Environment during calibration		
	Beginning	Finished
Temp. (°C)	24	24
REL.Humid. (%)	51	55
AC Supply (Volt)	223	223

Position :	Ref. Std. ID No.:
1	19-16RTD-01
2	19-16RTD-02
3	19-16RTD-03
4	23-16RTD-04
5	19-16RTD-05
6	19-14RTD-06
7	21-16RTD-07
8	19-16RTD-08
9 (ref.)	19-16RTD-09

Probe Installation Details :

a = 10 cm
b = 10 cm
c = 10 cm

Dimension of Chamber :

D = 0.50 m
W = 0.75 m
H = 1.2 m
Capacity = 0.45 m³



Equipment : Incubator
Condition As-Received : Used Item
Reference : 2412-0004OC-1
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Close

Cert. No.: 24TM1398
Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor k
35.0	35.0	35.0	0.048	0.40	0.46	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	
35.0	34.888	34.840	35.116	35.141	34.750	34.896	34.921	35.054	34.768	0.30

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor **k**, providing a level of confidence of approximately 95 %.

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert. No.: 24TM667
Page : 1 of 3

Equipment : Hot Air Oven
Manufacturer : Binder
Model : ED 240/E2
Serial No. : 00-15533
ID No. : BKK_ML0013

Submitted by : ALS Laboratory Group (Thailand) Co.,Ltd.
104 Phatthanakan 40, Phatthanakan Rd.,
Khwaeng Phatthanakan, Khet Suan Luang,
Bangkok 10250 Thailand

Location : Media Preparation Room

Received Order : 23 April 2024
Calibration Date : 23 April 2024
Ambient Temperature : $(26 \pm 10) ^\circ\text{C}$
Relative Humidity : $(50 \pm 30) \%$

Calibrated by : Tawatchai Pama

Approved by :

- () Ponpan Paipim
(✓) Suwit Imjai
() Kunchit Promprat

Issue Date : 26 April 2024

REVIEW BY Sithichok T.
APPROVED BY [Signature]
NEXT CAL DATE 23/10/25

[Signature]
Approved Signatory

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Equipment : Hot Air Oven
Condition As-Received : Used Item
Reference : 2404-0439OC-8
Procedure Used :-

Cert. No.: 24TM667
Page : 2 of 3

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Thermocouple Type T.

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Data Acquisition	MY49001451	24LM44	TPA	17 Mar 2025
2. This certificate is valid only to the item calibrated on date and place of calibration.				
3. This certification is traceable to the International System of Unit.				


Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Close

Environment during calibration		
	Beginning	Finished
Temp. (°C)	24	23
REL.Humid. (%)	65	65
AC Supply (Volt)	223	222

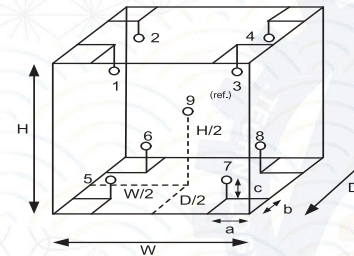

Probe Installation Details :

a = 10 cm
b = 10 cm
c = 10 cm

Dimension of Chamber :

D = 0.50 m
W = 0.80 m
H = 0.60 m
Capacity = 0.24 m³

Position :	Ref. Std. ID No.:
1	24-19TC-01
2	24-19TC-02
3	24-19TC-03
4	24-19TC-04
5	24-19TC-05
6	24-19TC-06
7	24-19TC-07
8	24-19TC-08
9 (ref.)	24-19TC-09



Equipment : Hot Air Oven
Condition As-Received : Used Item
Reference : 2404-0439OC-8
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Close

Cert. No.: 24TM667
Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor <i>k</i>
180	180	180	0.64	2.7	3.7	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	
180	181.009	181.511	180.922	181.359	181.217	183.659	181.664	181.986	181.474	1.5

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-oOo-



Internal Calibration

REVIEW BY	Somsak J.
APPROVED BY	Kanitta H.
NEXT CAL. DATE	22 May 25

Equipment : Turbidity Meter Manufacture : HACH
ID No. : SGK_FS0045 Model : 2100Q/QIS
Calibrate Date : May 22, 2024 Serial No. : 19010C073443

Calibration Point	1st (NTU)	2nd (NTU)	3rd (NTU)	AVG (NTU)	Specifications	Evaluate
Standard 20 NTU	19.8	20.1	20.0	20.0	19 to 21 NTU	Pass
Standard 100 NTU	99.5	99.3	100	99.6	95 to 105 NTU	Pass
Standard 800 NTU	798	799	798	798.3	760 to 840 NTU	Pass

Calibrated by Somsak J.
(Scientist 2)

Approved by : Kanitta H.
(Section Head)

HACH COMPANY



An ISO 9001 Certified Company

P.O.Box 389
Loveland, CO 80539
(970) 669-3050

Certificate of Analysis

This is a Component of 2659405-TH lot A4025

Page 1

COMMODITY: STABLCAL|sup|TS|sup0 FORMAZIN STANDARD 20 NTU
COMMODITY NUMBER: 2684801 MANUFACTURE DATE: DATE OF ANALYSIS:
LOT NUMBER: A4025 2/6/2024 2/12/2024

TEST	SPECIFICATIONS	RESULTS
Turbidity	19 to 21 NTU	20.5 NTU

The expiration date is May 2025

Formazin and StablCal® solutions provided by Hach are not NIST traceable because the NIS does not carry turbidity standards. However, the use of Formazin and StablCal® as use in Hach method 8195 are accepted by the EPA as a primary standard to be used in the calibration of turbidity instruments.

Certified by

Scott Als

Scott Als
Analytical Services Chemist

HACH COMPANY



An ISO 9001 Certified Company

P.O.Box 389
Loveland, CO 80539
(970) 669-3050

Certificate of Analysis

This is a Component of 2659405-TH lot A4025

Page 1

COMMODITY: STABLCAL|sup|TS|sup0 FORMAZIN STANDARD 100 NTU
COMMODITY NUMBER: 2684901 MANUFACTURE DATE: DATE OF ANALYSIS:
LOT NUMBER: A4018 1/29/2024 2/2/2024

TEST	SPECIFICATIONS	RESULTS
Turbidity	95 to 105 NTU	102.0 NTU

The expiration date is Apr 2025

Formazin and StablCal® solutions provided by Hach are not NIST traceable because the NIS does not carry turbidity standards. However, the use of Formazin and StablCal® as use in Hach method 8195 are accepted by the EPA as a primary standard to be used in the calibration of turbidity instruments.

Certified by

Scott Als

Scott Als
Analytical Services Chemist

HACH COMPANY



An ISO 9001 Certified Company

P.O.Box 389
Loveland, CO 80539
(970) 669-3050

Certificate of Analysis

This is a Component of 2659405-TH lot A4025

Page 1

COMMODITY: STABLCAL|sup|TS|sup0 FORMAZIN STANDARD 800 NTU
COMMODITY NUMBER: 2660501 MANUFACTURE DATE: DATE OF ANALYSIS:
LOT NUMBER: A4017 2/1/2024 2/9/2024


TEST	SPECIFICATIONS	RESULTS
Turbidity	760 to 840 NTU	818.0 NTU

The expiration date is May 2025

Formazin and StablCal® solutions provided by Hach are not NIST traceable because the NIS does not carry turbidity standards. However, the use of Formazin and StablCal® as use in Hach method 8195 are accepted by the EPA as a primary standard to be used in the calibration of turbidity instruments.

Certified by

Scott Als
Analytical Services Chemist

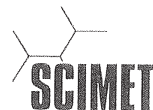

SCIMET Co., Ltd.
1194 Soi Wachirathamsathit 57, Bangchak,
Phrakhanong, Bangkok 10260 Thailand
Email: scimet2022@gmail.com, Tel: 02 460 9239
https://www.scimet.co.th



Certificate No. C07240192

Calibration Certificate

Equipment: SPECTROPHOTOMETER
Model: DR 3900
Serial No.(or ID): 2403637
Manufacturer: HACH
Condition: New

Job No.: KSMT2403543
Received Date: 25 December 2024
Issued Date: 25 December 2024
Page: 1 of 3

Customer

ALS Laboratory Group (Thailand) Co., Ltd
114/1 Moo 8, Karnchanawanich Rd.T. Ban Phru, A Hatyai, Songkhla 90250

Calibration Place

Hach (Thailand) Limited.
Branch 00001, Building D Room No. D3 11, 3rd Floor, No. 735/4, Srinakarin Road,
Pattanakarn, Suanluang, Bangkok 10250 Thailand.

Calibration Date

25 December 2024

Environment Condition

Temperature: 22.2 °C ± 0.5 °C
Humidity: 59.2 %RH ± 2.3 %RH

The Method used

In-house method, WI07, based on ASTM E 275-08 and
ASTM E 387-04

Traceability

This certificate is traceable to the CRM maintained by National Institute of Standards and Technology (NIST) through Starna Scientific Limited.

The standard for Wavelength Certificate No. 108691 and 108692

The standard for Photometric Certificate No. 109010, 114655

This certificate is issued the units of measurement according to the International System of Units (SI). It provides traceability of measurement to international or national standard or other recognized national standard laboratories.

The measurement uncertainty stated is the expanded uncertainty which is obtained from the standard uncertainty multiplied by the coverage factor ($k=2$) to provide a level of confidence of approximately 95%. It is determined in accordance with the Guide to Expression of Uncertainty in Measurement (GUM).

These results may be affected by deviations from specified conditions. The results relate only to the items tested, calibrated or sampled. The report shall not be reproduced except in full without approval of SCIMET Co., Ltd.

(Mr. Siwapan Srijan)

Person in charge



(Mr. Thalerngkeat Pongngam)

Authorized signatory



Certificate No.: C07240192 Page 2 of 3

Condition of reference standards Instruments / CRM:

Instruments	Set No.	Certificate No.	Due date
Holmium Oxide Glass Reference	121512	108691	25-Jan-25
Didymium Oxide Glass Reference	119722	108692	25-Jan-25
Neutral Density Filter Reference	12276	109010 , 114655	2-Feb-25

Calibration Results:

Without Adjustment

Wavelength Accuracy (nm), The spectral bandwidth of Std at 5 nm and UUC at 5 nm

Standard Wavelength (nm)	Unit Under Calibration (nm)	Correction (nm)	Uncertainty of Measurement (± nm)
333.67	333	0.67	0.59
361.02	361	0.02	0.59
417.80	417	0.80	0.59
441.29	441	0.29	0.59
479.88	480	-0.12	0.59
513.75	513	0.75	0.59
528.59	528	0.59	0.59
537.75	537	0.75	0.59
585.56	585	0.56	0.59
641.95	642	-0.05	0.59
684.70	685	-0.30	0.59
747.61	748	-0.39	0.59
807.04	807	0.04	0.59
879.68	880	-0.32	0.59



Certificate No.: C07240192 Page 3 of 3

Calibration Results:

Without Adjustment

Photometric Accuracy (Absorbance)

Wavelength	Standard absorbance (Abs)	Unit Under Calibration (Abs)	Correction (Abs)	Uncertainty of Measurement(± Abs)
420 nm	0.0000	0.000	0.0000	0.0045
	0.2373	0.234	0.0033	0.0045
	0.5617	0.561	0.0007	0.0045
	0.7392	0.738	0.0012	0.0045
	1.0550	1.056	-0.0010	0.0045
440 nm	0.0000	0.000	0.0000	0.0045
	0.2335	0.231	0.0025	0.0045
	0.5513	0.550	0.0013	0.0045
	0.7230	0.722	0.0010	0.0045
	1.0324	1.032	0.0004	0.0045
465 nm	0.0000	0.000	0.0000	0.0045
	0.2126	0.210	0.0026	0.0045
	0.5036	0.506	-0.0024	0.0045
	0.6735	0.676	-0.0025	0.0045
	0.9615	0.965	-0.0035	0.0045
546.1 nm	0.0000	0.000	0.0000	0.0045
	0.2201	0.217	0.0031	0.0045
	0.5176	0.520	-0.0024	0.0045
	0.6930	0.694	-0.0010	0.0045
	0.9908	0.994	-0.0032	0.0045
590 nm	0.0000	0.000	0.0000	0.0045
	0.2443	0.241	0.0033	0.0045
	0.5530	0.554	-0.0010	0.0045
	0.7196	0.719	0.0006	0.0045
	1.0301	1.031	-0.0009	0.0045
635 nm	0.0000	0.000	0.0000	0.0045
	0.2646	0.261	0.0036	0.0045
	0.5370	0.538	-0.0010	0.0045
	0.6862	0.687	-0.0008	0.0045
	0.9822	0.984	-0.0018	0.0045

The End of Certificate

ภาคผนวก จ

สำเนาหนังสืออนุญาตขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
จากกรมโรงงานอุตสาหกรรม



ที่อก ๐๓๑๐(๑)/ ๑๖๑๖๘

กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๕๐๐

๒๐ พฤศจิกายน ๒๕๖๖

เรื่อง ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

เรียน กรรมการผู้จัดการ บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๔ สิงหาคม ๒๕๖๖

- สิ่งที่ส่งมาด้วย ๑. รายชื่อผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๑ แผ่น
๒. รายชื่อเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๕ แผ่น
๓. ขอบข่ายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๓๑ แผ่น

ตามหนังสือที่อ้างถึง บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด ขอต่ออายุหนังสือ
รับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ว-๒๐๔ สถานที่ตั้งเลขที่ ๑๐๔ ซอยพัฒนาการ ๔๐
ถนนพัฒนาการ แขวงพัฒนาการ เขตสวนหลวง กรุงเทพมหานคร ต่อกรมโรงงานอุตสาหกรรม นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด
ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน โดยมีองค์ประกอบดังนี้

- ก. ผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๖ ราย ตามสิ่งที่ส่งมาด้วย ๑
ข. เจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ ๑๘๑ ราย ตามสิ่งที่ส่งมาด้วย ๒
ค. ขอบข่ายสารมลพิษที่ได้รับขึ้นทะเบียนให้วิเคราะห์ในน้ำเสีย น้ำใต้ดิน อากาศเสีย สิ่งปฏิกูล
หรือวัสดุที่ไม่ใช้แล้ว และดิน ตามสิ่งที่ส่งมาด้วย ๓

หนังสือฉบับนี้จะหมดอายุในวันที่ ๒ กันยายน ๒๕๖๙ หากประสงค์จะต่ออายุหนังสือ
รับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ให้ยื่นคำขอต่ออายุพร้อมเอกสารประกอบคำขอต่อ
กรมโรงงานอุตสาหกรรม ภายใน ๓๐ วัน ก่อนวันสิ้นอายุของหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
ทั้งนี้ สามารถยื่นคำขอผ่านระบบอิเล็กทรอนิกส์ได้ที่หน้าเว็บไซต์กรมโรงงานอุตสาหกรรม

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

(นายสิระ จันทระดิษฐ์)

นักวิทยาศาสตร์เชี่ยวชาญ รักษาราชการแทน
ผู้อำนวยการกองวิจัยและเตือนภัยมลพิษโรงงาน
ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม

กองวิจัยและเตือนภัยมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕

โทรสาร ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๔๔

ไปรษณีย์อิเล็กทรอนิกส์ saraban@diw.mail.go.th



เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด เลขทะเบียน ว-๒๐๔
ที่ ออก ๐๓๑๐(๑)/ ๑ ๖ ๑ ๖ ๘ ลงวันที่ ๒๐ พฤศจิกายน ๒๕๖๖

ก. ผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์ จำนวน ๖ ราย

- | | |
|--------------------------------|----------------------------|
| ๑) นางสาวยุพาพร จันทรเปล่ง | ทะเบียนเลขที่ ว-๒๐๔-ค-๐๐๐๑ |
| ๒) นางสาวชัชชัย โกมารกุล ณ นคร | ทะเบียนเลขที่ ว-๒๐๔-ค-๐๐๐๒ |
| ๓) นายศรายุทธ จิตรานนท์ | ทะเบียนเลขที่ ว-๒๐๔-ค-๐๐๐๓ |
| ๔) นางสาวกนกกร เอนก | ทะเบียนเลขที่ ว-๒๐๔-ค-๐๐๐๔ |
| ๕) นายสุริยา สอนแก้ว | ทะเบียนเลขที่ ว-๒๐๔-ค-๐๐๐๕ |
| ๖) นายวิชาญ ชุณหะวัณ | ทะเบียนเลขที่ ว-๒๐๔-ค-๐๐๐๖ |

3/กข

เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด

เลขทะเบียน ว-๒๐๔

ที่ ออก ๐๓๑๐(๑)/ ๑๖๑๖๘

ลงวันที่ ๒๐ พฤศจิกายน ๒๕๖๖

ข. เจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๑๘๑ ราย

๑) นายกาจบัณฑิต กิตติคุณวิชัย	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๐๑
๒) นายภัทรพล สว่างใจธรรม	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๐๒
๓) นายณราธิป เทือกชัยคำ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๐๓
๔) นายศิริโชค พงษ์ประสม	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๐๔
๕) นายณัฐวุฒิ ดั่งแพง	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๐๕
๖) นางสาวจินดา ไชจุลธรรม	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๐๖
๗) นางสาวสาวิตรี น้อยเสงี่ยม	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๐๗
๘) นางสาวชนัญญาญจน์ อัมมขม	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๐๘
๙) นางสาวนรินทร์ สายเส็ง	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๐๙
๑๐) นางสาวนันทาศิ สมบูรณ์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๑๐
๑๑) นางสาวศรณิยา เฉลิมธารงค์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๑๑
๑๒) นางสาวธัญญธร มงคลจิรวุฒิ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๑๒
๑๓) นางสาวศิริลักษณ์ บุญนาค	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๑๓
๑๔) นายณพนงค์ จันทุพันธุ์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๑๔
๑๕) นายนรเศรษฐ์ โกมลาลัย	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๑๕
๑๖) นายธันวา จริยา	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๑๖
๑๗) นางสาวเกศรินทร์ แก้วมัน	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๑๗
๑๘) นางสาวสุวิมล ชัยเรืองวุฒิ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๑๘
๑๙) นางสาวสุชาดา ธรรมถาวร	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๑๙
๒๐) นางสาวเปมิกา ชัยเดชธนกุล	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๒๐
๒๑) นางสาวศศิธร หมูสวัสดิ์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๒๑
๒๒) นางสาวเสาวลักษณ์ ภูณภาพอำพร	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๒๒
๒๓) นายอภิสิทธิ์ สิงหา	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๒๓
๒๔) นายศักดิ์สิทธิ์ ไพศาลพิสุทธิ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๒๔
๒๕) ว่าที่ร้อยตรีหญิง พรรณิภา ขำเจริญ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๒๕
๒๖) นางจิตดา คำภูแก้ว	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๒๖
๒๗) นางสาวอรรพรรณ รักยง	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๒๗
๒๘) นางสาวนพรัตน์ แยมกรานต์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๒๘
๒๙) นายจุลเดช วารินทร์	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๒๙
๓๐) นางสาวดาญรัตน์ ร้องคำ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๓๐
๓๑) นายพรมมี ศรีปัตเนตร	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๓๑
๓๒) นายอุทิศ อุ่นสิม	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๓๒
๓๓) ว่าที่ร้อยตรี เฉลิมเกียรติ อมรศรีเสริม	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๓๓
๓๔) นางสาววริยา สร้างนา	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๓๔
๓๕) นายอนุพงศ์ รัตนศรีประเสริฐ	ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๓๕

๓๖) นางสาวจุฑารัตน์...

๓๖) นางสาวจุฑารัตน์ โอนสันเทียะ
๓๗) นางสาวจารุวรรณ พิมพ์ภักฤติยา
๓๘) นางสาวปรารถน์ทิพย์ กิจไพศาลศักดิ์
๓๙) นางสาวเดือนใจ ทางกลาง
๔๐) นางสาวจิราพร ศิริเวช
๔๑) นายวรากร ผูกרך
๔๒) นายทง วิริยะสทกิจ
๔๓) นายธนิศ เจนจบ
๔๔) นายคณิศร ขำเพชร
๔๕) นายภูวิช พรหมสะอาด
๔๖) นายธนเดช โกศาพิพัฒน์
๔๗) นายชวฤทธิ์ วงษ์จันทร์
๔๘) นายอาทิตย์ ศรีเสน
๔๙) นายเจษดินทร์ คงศักดิ์ไทย
๕๐) นายจรัส บุญยั้ง
๕๑) นายธนาณัติ เอนก
๕๒) นายอภิวัฒน์ ทุมหนู
๕๓) นางสาวสุภาขวัญ มาก
๕๔) นางสาวหัตถพร ขวาลสมบุรณ์
๕๕) นางสาวธิดิมา บุญเพ็ง
๕๖) นางสาวภานุมาศ นามวัฒน์
๕๗) นางสาวอุไรรัตน์ ทังสร้างแป้น
๕๘) นายธีรวัฒน์ ปวงสุข
๕๙) นายอิทธิพล ยะโส
๖๐) นายประพจน์ วรรณชูชัย
๖๑) นายชยธร พวงทิพย์
๖๒) นางสาวกนกวรรณ จันทบาล
๖๓) นายสิทธิโชค ธงเงิน
๖๔) นางศิวาวรรณ ใจบุญ
๖๕) นางสาวพรรณธิดา พุ่มคง
๖๖) นายนวกัทร ศรีวิริยะ
๖๗) นายสุริษา ทองอ่อน
๖๘) นายวิญญู บุญตะนัย
๖๙) นายสมบุรณ์ บุตรจันทร์
๗๐) นายวิรัตน์ ไชยนระรา
๗๑) นายณฤเบศน์ เพิ่มพูน
๗๒) นายจิรณัฐ ขาวละออ
๗๓) นายอัสนี นามบุรี
๗๔) นายอัศวเรศ จ่อสา

ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๓๖
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๓๗
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๓๘
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๓๙
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๔๐
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๔๑
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๔๒
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๔๓
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๔๔
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๔๕
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๔๖
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๔๗
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๔๘
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๔๙
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๕๐
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๕๑
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๕๒
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๕๓
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๕๔
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๕๕
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๕๖
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๕๗
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๕๘
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๕๙
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๐
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๑
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๒
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๓
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๔
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๕
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๖
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๗
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๘
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๙
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๗๐
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๗๑
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๗๒
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๗๓
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๗๔

๗๕) นายประเสริฐ สุระขันธุ์
๗๖) นายอนุกุล จันทน์เนียม
๗๗) นายพิรพงษ์ ทองคุณปรีดา
๗๘) นายณฤพล ทองนุช
๗๙) นายอนุวัฒน์ ม่วงแพ
๘๐) นายเจตศราวุฒิ ปัตตะมะ
๘๑) นายกฤษณะ สายวรรณ
๘๒) นายพิชัย บุญยงค์
๘๓) นายภาณุพงศ์ โหมวงศ์
๘๔) นายสามารถ คุ่มปลี
๘๕) นายสัญญาชัย โกศรีนาม
๘๖) นายณัฐวุฒิ ศรีประเสริฐ
๘๗) นายชวัลรัช นาคพนม
๘๘) นายพงศธร ชัยทิพย์
๘๙) นายสิทธิโชค ทาสีดา
๙๐) นายธนากร อินสุตา
๙๑) นางสาววรรณิษา ชาตวันชัย
๙๒) นางสาวพิมพ์ตะวัน มินากุล
๙๓) นางสาวเพชรรัตน์ สิงห์สมบุญ
๙๔) นางสาวชญาณิน พรหมจันทร์
๙๕) นายกীরติ ทวีราช
๙๖) นายจักริน หมั่นวิชา
๙๗) นายฉัตรชัย สุขเปี้ย
๙๘) นายณรรนท ต๊ะทองคำ
๙๙) นายศุภผล สนนอก
๑๐๐) นายทักษ์ดนัย อุบลศรี
๑๐๑) นายธนากร นามะกฤษณา
๑๐๒) นายธิตพงศ์ บัวแดง
๑๐๓) นายณนทชัย อุปถัมภ์
๑๐๔) นายณัฐพล คุณสุทธิ
๑๐๕) นายณัฏฐวัฒน์ สาริน
๑๐๖) นายปิยะนัฐ พลมะศรี
๑๐๗) นายพงศ์สิริ โสมเขียว
๑๐๘) นายพีรพัฒน์ กำคำ
๑๐๙) นายภาณุพงศ์ มานิตย์
๑๑๐) นายมงคล ผลาทิพย์
๑๑๑) นายสิรินนท์ ทองอ้น
๑๑๒) นายอเนชา ทันสมัย
๑๑๓) นายอดิศักดิ์ ผมไผ

ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๗๕
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๗๖
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๗๗
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๗๘
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๗๙
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๘๐
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๘๑
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๘๒
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๘๓
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๘๔
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๘๕
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๘๖
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๘๗
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๘๘
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๘๙
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๙๐
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๙๑
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๙๒
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๙๓
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๙๔
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๙๕
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๙๖
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๙๗
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๙๘
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๙๙
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๐๐
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๐๑
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๐๒
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๐๓
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๐๔
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๐๕
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๐๖
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๐๗
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๐๘
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๐๙
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๑๐
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๑๑
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๑๒
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๑๓

วิมล

๑๑๔) นายอนันตชัย...

๑๑๔) นายอนันตชัย วิสม
๑๑๕) นายวรวิธ ตีนัก
๑๑๖) นายแสงตะวัน นตะสัด
๑๑๗) นายยุทธพงศ์ รัตนะ
๑๑๘) นายชัยวัฒน์ ไชยชนะนิจ
๑๑๙) นายวิศรุต ศรีธรรมมา
๑๒๐) นายณนทกร เผือกผ่อง
๑๒๑) นายกำชัย สุทธะ
๑๒๒) นางสาวณัฐภรณ์ บุญตะนัย
๑๒๓) นางสาวพัชรินทร์ แสนสร้อย
๑๒๔) นายไพรวลัย เปี่ยมพิมาย
๑๒๕) นางสาวศุภมาศ ทองมาก
๑๒๖) นางสาวลลิตา จิตรสว่าง
๑๒๗) นางสาวไมพร เสิกภูเขียว
๑๒๘) นางสาวกฤติมาพร คำมีแก่น
๑๒๙) นางสาวสกุลรัตน์ ภาคภูมิ
๑๓๐) นางสาวไพรินทร์ ศรีรูปี
๑๓๑) นางสาวทิพนเตร ผุยปัญญา
๑๓๒) นางสาวสาธิตา ปานทอง
๑๓๓) นางสาวอริสา ทองนวล
๑๓๔) นางสาวอรยา คำคล้อง
๑๓๕) นางสาวชุตาภรณ์ สุนทรสนาน
๑๓๖) นางสาวอัญชลี คำจันทร์
๑๓๗) นายบุญฤทธิ์ เอี่ยมเทศ
๑๓๘) นางสาวศุภรดา บันมยุรา
๑๓๙) นางสาวพาฤดี คุณนาน
๑๔๐) นางสาวจิราเจต ฟองดา
๑๔๑) นางสาวอารยา มีชัย
๑๔๒) นางสาววิษุตา นาคผจญ
๑๔๓) นางสาวนันทิยา จันทะสุน
๑๔๔) นายกิตติพงศ์ แซ่ลี
๑๔๕) นายอนวัติ ภูถวิล
๑๔๖) นายธีรพล แสงทอง
๑๔๗) นายศักดิ์พิพัฒน์ บุญมัน
๑๔๘) นายฐิติวัลด์ เอมอุไร
๑๔๙) นายชัยณรงค์ ศรีบุรินทร์
๑๕๐) นางสาวอัจฉราวรรณ สอนสนอง
๑๕๑) นางสาวณัฐพร สิงหา
๑๕๒) นายกัมเรศ แหยมโต

ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๑๔
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๑๕
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๑๖
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๑๗
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๑๘
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๑๙
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๒๐
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๒๑
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๒๒
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๒๓
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๒๔
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๒๕
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๒๖
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๒๗
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๒๘
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๒๙
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๓๐
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๓๑
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๓๒
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๓๓
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๓๔
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๓๕
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๓๖
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๓๗
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๓๘
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๓๙
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๔๐
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๔๑
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๔๒
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๔๓
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๔๔
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๔๕
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๔๖
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๔๗
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๔๘
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๔๙
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๕๐
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๕๑
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๕๒

๑๕๓) นางสาวอุบล เคิกศิริ
๑๕๔) นางสาวมโนรัตน์ ทองบุตร
๑๕๕) นายภาคภูมิ แทนไทย
๑๕๖) นางสาวสุภาณัฐ เมล์พวง
๑๕๗) นางสาวพรทิศา สาดาชนม์
๑๕๘) นายเอกวิทย์ วันทะนา
๑๕๙) นายไตรมณฑล ทิพย์วรรณ
๑๖๐) นายจิรเมธ ประเสริฐสิริพงศ์
๑๖๑) นายจิรายุส เกษมสุข
๑๖๒) นายจิรศักดิ์ ศรีวิชัย
๑๖๓) นายณัฐกฤษณ์ สะพานแก้ว
๑๖๔) นายบุญศักดิ์ ปะที
๑๖๕) นายปณณวิทย์ เสมอทรัพย์
๑๖๖) นายพิษณุพงษ์ ไชยา
๑๖๗) นายภัทรพงษ์ มณฑาทอง
๑๖๘) นายวสันต์ ตรีนกุล
๑๖๙) นายภาณุเดช เพชรอุด
๑๗๐) นายอนุกุล วิละแสง
๑๗๑) นายภัทรพงษ์ มีสุข
๑๗๒) นางสาวนุชรี ลีละทีป
๑๗๓) นางสาวสุภาวดี โกศรีนาม
๑๗๔) นางสาวอรณิข เทียนคำ
๑๗๕) นางสาวพรเพ็ญ ชอบสอน
๑๗๖) นางสาววันวิสา ขอนพิกุล
๑๗๗) นางสาวอรรณณ เถาว์ทอง
๑๗๘) นางสาวอัยยลิณ เมอร์วิณณ์
๑๗๙) นางสาววิสรา ค่วยครอง
๑๘๐) นายวุฒิกร ศิริวรรณ
๑๘๑) นางสาวจรรวรรณ กระจำงพันธุ์

ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๕๓
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๕๔
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๕๕
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๕๖
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๕๗
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๕๘
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๕๙
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๐
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๑
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๒
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๓
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๔
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๕
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๖
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๗
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๘
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๙
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๗๐
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๗๑
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๗๒
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๗๓
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๗๔
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๗๕
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๗๖
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๗๗
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๗๘
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๗๙
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๐
ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๑

วิมล

เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท เอแอลเอส แลบบอราทอรี กรู๊ป (ประเทศไทย) จำกัด เลขทะเบียน ว-๒๐๔

ที่ อก ๐๓๑๐(๑)/ ๑๖๑๖๘ ลงวันที่ ๒๐ พฤศจิกายน ๒๕๖๖

ค. ขอบข่ายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๓๗๔ รายการ

น้ำเสีย จำนวน 60 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Aldicarb	High-Performance Liquid Chromatographic Method ^[4]
2	Aldicarb Sulfone	High-Performance Liquid Chromatographic Method ^[4]
3	Aldicarb Sulfoxide	High-Performance Liquid Chromatographic Method ^[4]
4	Aldrin	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
5	Arsenic	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
6	Barium	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
7	α -BHC	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
8	β -BHC	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
9	δ -BHC	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
10	γ -BHC	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
11	Biochemical Oxygen Demand	1) 5-Day BOD Test, Azide Modification Method ^[4] 2) 5-Day BOD Test, Membrane Electrode Method ^[4]
12	Carbaryl	High-Performance Liquid Chromatographic Method ^[4]
13	Carbofuran	High-Performance Liquid Chromatographic Method ^[4]
14	Cadmium	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
15	Chemical Oxygen Demand	1) Closed Reflux, Colorimetric Method ^[4] 2) Closed Reflux, Titrimetric Method ^[4]
16	Chlordane	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
17	Chromium	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
18	Color	ADMI Weighted-Ordinate Spectrophotometric Method ^[4]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
19	Copper	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
20	Cyanide	Distillation, Colorimetric Method ^[4]
21	2,4'-DDD	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
22	4,4'-DDD	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
23	2,4'-DDE	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
24	4,4'-DDE	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
25	2,4'-DDT	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
26	4,4'-DDT	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
27	Dieldrin	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
28	Endosulfan Sulfate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
29	Endosulfan I	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
30	Endosulfan II	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
31	Endrin	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
32	Endrin Aldehyde	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
33	Formaldehyde	Distillation, Colorimetric Method ^[3]
34	Free Chlorine	1) DPD Ferrous Titrimetric Method ^[4] 2) DPD Colorimetric Method ^[4]
35	Heptachlor	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
36	Heptachlor Epoxide	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
37	Hexavalent Chromium	Colorimetric Method ^[4]
38	3-Hydroxycarbofuran	High-Performance Liquid Chromatographic Method ^[4]
39	Lead	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
40	Manganese	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
41	Mercury	1) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass spectrometric Method ^[4]
42	Methiocarb	High-Performance Liquid Chromatographic Method ^[4]
43	Methoxychlor	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
44	Methomyl	High-Performance Liquid Chromatographic Method ^[4]
45	Nickel	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
46	Oil & Grease	1) Liquid-Liquid, Partition-Gravimetric Method ^[4] 2) Soxhlet Extraction Method ^[4]
47	Oxamyl	High-Performance Liquid Chromatographic Method ^[4]
48	Propoxur	High-Performance Liquid Chromatographic Method ^[4]
49	pH	Electrometric Method ^[4]
50	Phenols	1) Distillation, Chloroform Extraction Method ^[4] 2) Distillation, Direct Photometric Method ^[4]
51	Selenium	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
52	Sulfide	Iodometric Method ^[4]
53	Temperature	Laboratory and Field Methods ^[4]
54	Total Dissolved Solids	Dried at 180 °C ^[4]
55	Total Kjeldahl Nitrogen	Semi-Micro Kjeldahl Method ^[4]
56	Total Phosphorous	Digestion, Colorimetric Method ^[4]
57	Total Suspended Solids	Dried from 103-105 °C ^[4]
58	Toxaphene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
59	Trivalent Chromium	1) Digestion, Inductively Coupled Plasma Method; Colorimetric Method; Calculation ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method; Colorimetric Method; Calculation ^[4]
60	Zinc	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[4]

น้ำใต้ดิน จำนวน 126 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Acenaphthene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
2	Acetone	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
3	Aldrin	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
4	Anthracene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
5	Antimony	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
6	Arsenic	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
7	Atrazine	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
8	Barium	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
9	Benz(a)anthracene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
10	Benzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
11	Benzo(b)fluoranthene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
12	Benzo(k)fluoranthene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
13	Benzoic Acid	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
14	Benzo(a)pyrene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
15	Benzo[g,h,i]perylene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
16	Beryllium	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
17	Bis(2-chloroethyl)ether	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
18	Bis(2-ethylhexyl)phthalate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
19	Bromodichloromethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
20	Bromoform	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
21	Butanol	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
22	Butyl benzyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
23	Cadmium	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
24	Carbazole	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
25	Carbon disulfide	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
26	Carbon tetrachloride	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
27	Chlordane	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
28	p-Chloroaniline	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
29	Chlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
30	Chlorodibromomethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
31	Chloroform	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
32	2-Chlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
33	Chromium	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
34	Chromium (III)	1) Digestion, Inductively Coupled Plasma Method; Colorimetric Method; Calculation ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method; Colorimetric Method; Calculation ^[4]
35	Chromium (VI)	Colorimetric Method ^[4]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
36	Chrysene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
37	Cyanide	Distillation, Colorimetric Method ^[4]
38	2,4-D	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
39	DDD	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
40	DDE	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
41	DDT	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
42	Dibenz(a,h)anthracene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
43	Di-n-Butyl Phthalate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
44	1,2-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
45	1,3-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
46	1,4-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
47	3,3-Dichlorobenzidine	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
48	1,1-Dichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
49	1,2-Dichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
50	1,1-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
51	cis-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
52	trans-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
53	2,4-Dichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
54	1,2-Dichloropropane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
55	1,3-Dichloropropane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
56	1,3-Dichloropropene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
57	Dieldrin	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
58	Diethyl Phthalate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
59	2,4-Dimethylphenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
60	2,4-Dinitrophenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
61	2,4-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
62	2,6-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
63	Di-n-octyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
64	Endosulfan	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
65	Endrin	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
66	Ethylbenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
67	Fluoranthene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
68	Fluorene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
69	Heptachlor	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
70	Heptachlor epoxide	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
71	Hexachlorobenzene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
72	Hexachloro-1,3-butadiene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
73	n-Hexane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
74	α-HCH	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
75	β-HCH	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
76	γ-HCH	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
77	Hexachlorocyclopentadiene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
78	Hexachloroethane	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
79	Indeno(1,2,3-cd)pyrene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
80	Isophorone	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
81	Lead	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
82	Manganese	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
83	Mercury	1) Digestion, Cold Vapor Atomic Absorption Spectrometric Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
84	Methanol	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
85	Methoxychlor	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
86	Methyl bromide	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
87	Methylene chloride	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
88	2-Methylphenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
89	2-Methylnaphthalene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
90	Methyl tert-butyl Ether	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
91	Naphthalene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
92	Nickel	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
93	Nitrobenzene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
94	N-Nitrosodiphenylamine	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
95	N-Nitrosodi-n-Propylamine	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
96	Polychlorinated Biphenyls - PCB 1016 - PCB 1221 - PCB 1232 - PCB 1242 - PCB 1248 - PCB 1254 - PCB 1260	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
97	Pentachlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
98	pH	Electrometric Method ^[4]
99	Phenanthrene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
100	Phenol	1) Distillation, Chloroform Extraction Method ^[4] 2) Distillation, Direct Photometric Method ^[4] 3) Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
101	Pyrene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
102	Selenium	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
103	Silver	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[4]
104	Styrene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
105	1,1,2,2-Tetrachloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
106	Tetrachloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
107	Toluene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[4]
108	Toxaphene	Liquid-Liquid Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[4]
109	TPH (C ₅ -C ₈)	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[14,25]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
110	TPH (C ₈ -C ₁₆)	Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[9,22]
111	TPH (C ₁₆ -C ₃₅)	Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[9,22]
112	1,2,4-Trichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
113	1,1,1-Trichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
114	1,1,2-Trichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
115	Trichloroethylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
116	2,4,5-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
117	2,4,6-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[4]
118	1,3,5-Trimethylbenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
119	Vanadium	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[4]
120	Vinyl acetate	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
121	Vinyl chloride	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
122	m-Xylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
123	o-Xylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
124	p-Xylene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
125	Xylene (Total)	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[4]
126	Zinc	1) Digestion, Inductively Coupled Plasma Method ^[4] 2) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[4]

อากาศเสีย (ปล่อยระบาย) จำนวน 28 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Antimony	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
2	Arsenic	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
3	Beryllium	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
4	Cadmium	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
5	Carbon Monoxide	1) Instrumental Analyzer Method ^[5] 2) Sampling Bag Non-Dispersive Infrared Method ^[5]
6	Chlorine	1) Absorption Sampling, Ion Chromatographic Method ^[5] 2) Isokinetic Sampling, Ion Chromatographic Method ^[5]
7	Chromium	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
8	Cobalt	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
9	Copper	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
10	Cresol	Adsorption Sampling, Gas Chromatographic Method ^[5]
11	Dioxins	Isokinetic Sampling ^[5]
12	Hydrogen Chloride	1) Absorption Sampling, Ion Chromatographic Method ^[5] 2) Isokinetic Sampling, Ion Chromatographic Method ^[5]
13	Hydrogen Fluoride	1) Absorption Sampling, Ion Chromatographic Method ^[5] 2) Isokinetic Sampling, Ion Chromatographic Method ^[5]
14	Hydrogen Sulfide	Absorption Sampling, Iodometric Method ^[5]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
15	Lead	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
16	Manganese	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
17	Mercury	1) Isokinetic Sampling, Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[5] 2) Isokinetic Sampling, Digestion, Cold-Vapor Atomic Fluorescence Spectrometric Method ^[5]
18	Nickel	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
19	Opacity	Ringelmann's Method ^[2]
20	Oxides of Nitrogen	1) Absorption Sampling, Phenoldisulfonic Acid Method ^[5] 2) Absorption Sampling, Alkaline Permanganate/Colorimetric Method ^[5] 3) Instrumental Analyzer Method ^[5]
21	Selenium	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
22	Sulfur Dioxide	1) Absorption Sampling, Barium-Thorin Titrimetric Method ^[5] 2) Instrumental Analyzer Method ^[5]
23	Sulfuric Acid	Isokinetic Sampling, Barium-Thorin Titrimetric Method ^[5]
24	Tellurium	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
25	Tin	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
26	Total Suspended Particulate	1) Isokinetic Sampling, Gravimetric Method ^[5] 2) Paired Train, Isokinetic Sampling, Gravimetric Method ^[5]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
27	Vanadium	1) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[5]
28	Xylene	Adsorption Sampling, Gas Chromatographic Method ^[5]

สิ่งปฏิกูลหรือวัสดุที่ไม่ใช้แล้ว จำนวน 35 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Aldrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26]
2	Antimony	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]
3	Arsenic	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]
4	Barium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
5	Beryllium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]
6	Cadmium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]
7	Chlordane	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26]
8	Chromium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]
9	Chromium (III)	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method; Waste Extraction, Colorimetric Method; Calculation Method ^[1,6,16,19] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method; Waste Extraction, Colorimetric Method; Calculation Method ^[1,6,17,19] 3) Digestion, Inductively Coupled Plasma Method; Alkaline Digestion, Colorimetric Method; Calculation Method ^[7,8,16,19] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method; Alkaline Digestion, Colorimetric Method; Calculation Method ^[7,8,17,19]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
10	Chromium (VI)	1) Waste Extraction, Colorimetric Method ^[1,6,19] 2) Alkaline Digestion, Colorimetric Method ^[8,19]
11	Cobalt	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]
12	Copper	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]
13	2,4-D	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26]
14	DDD	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26]
15	DDE	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26]
16	DDT	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
17	Dieldrin	2) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26] 1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
18	Endrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
19	Heptachlor	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
20	Lead	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]
21	Lindane	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
22	Mercury	1) Waste Extraction, Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[1,6,20] 2) Waste Extraction, Digestion, Cold-Vapor Atomic Fluorescence Spectrometric Method ^[1,6,30] 3) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[20] 4) Digestion, Cold-Vapor Atomic Fluorescence Spectrometric Method ^[30] 5) Thermal Decomposition Amalgamation and Atomic Absorption Spectrometric Method ^[21]
23	Methoxychlor	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic / Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic / Mass Spectrometric Method ^[11,26]
24	Mirex	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic / Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic / Mass Spectrometric Method ^[11,26]
25	Molybdenum	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma / Mass Spectrometric Method ^[7,17]
26	Nickel	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma / Mass Spectrometric Method ^[7,17]
27	Polychlorinated biphenyls (PCBs) - Aroclor 1016 - Aroclor 1221 - Aroclor 1232 - Aroclor 1242 - Aroclor 1248 - Aroclor 1254 - Aroclor 1260	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic Method ^[11,26]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
	- 2-Chlorobiphenyl - 2,3-Dichlorobiphenyl - 2,2',5-Trichlorobiphenyl - 2,4',5-Trichlorobiphenyl - 2,2',3,5'-Tetrachlorobiphenyl - 2,2',5,5'-Tetrachlorobiphenyl - 2,3',4,4'-Tetrachlorobiphenyl - 2,2',3,4,5'-Pentachlorobiphenyl - 2,2',4,5,5'-Pentachlorobiphenyl - 2,3,3',4',6-Pentachlorobiphenyl - 2,2',3,4,4',5'-Hexachlorobiphenyl - 2,2',3,4,5,5'-Hexachlorobiphenyl - 2,2',3,5,5',6-Hexachlorobiphenyl - 2,2',4,4',5,5'-Hexachlorobiphenyl - 2,2',3,3',4,4',5-Heptachlorobiphenyl - 2,2',3,4,4',5,5'-Heptachlorobiphenyl - 2,2',3,4,4',5',6-Heptachlorobiphenyl - 2,2',3,4',5,5',6-Heptachlorobiphenyl - 2,2',3,3',4,4',5,5',6-Nonachlorobiphenyl	
28	Pentachlorophenol	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26]
29	pH	Electrometric Method ^[23,24]
30	Selenium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
31	Silver	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]
32	Thallium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]
33	Toxaphene	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^[1,9,26] 2) Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,26] 3) Automated Soxhlet Extraction, Gas Chromatographic/Mass Spectrometric Method ^[11,26]
34	Vanadium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]
35	Zinc	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[1,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1,6,17] 3) Digestion, Inductively Coupled Plasma Method ^[7,16] 4) Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[7,17]

ดิน จำนวน 125 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Acenaphthene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
2	Acetone	1) Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25] 2) Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method ^[13]
3	Aldrin	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
4	Anthracene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
5	Antimony	1) Digestion, Inductively Coupled Plasma Method ^[7,16] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]
6	Arsenic	1) Digestion, Inductively Coupled Plasma Method ^[7,16] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]
7	Atrazine	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
8	Barium	1) Digestion, Inductively Coupled Plasma Method ^[7,16] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]
9	Benz(a)anthracene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
10	Benzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
11	Benzo(b)fluoranthene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
12	Benzo(k)fluoranthene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
13	Benzoic acid	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
14	Benzo(a)pyrene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
15	Benzo(g,h,i)perylene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
16	Beryllium	1) Digestion, Inductively Coupled Plasma Method ^[7,16] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]
17	Bis(2-chloroethyl)ether	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
18	Bis(2-ethylhexyl)phthalate	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
19	Bromodichloromethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
20	Bromoform	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
21	Butanol	Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
22	Butyl Benzyl Phthalate	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
23	Cadmium	1) Digestion, Inductively Coupled Plasma Method ^[7,16] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]
24	Carbazole	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
25	Carbon Disulfide	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
26	Carbon tetrachloride	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
27	Chlordane	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
28	p-Chloroaniline	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
29	Chlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
30	Chlorodibromomethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
31	Chloroform	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
32	2-Chlorophenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
33	Chromium	1) Digestion, Inductively Coupled Plasma Method ^[7,16] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]
34	Chromium (III)	1) Digestion, Inductively Coupled Plasma Method; Alkaline Digestion, Colorimetric Method; Calculation Method ^[7,8,16,19] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method; Alkaline Digestion, Colorimetric Method; Calculation Method ^[7,8,17,19]
35	Chromium (VI)	Alkaline Digestion, Colorimetric Method ^[8,19]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
36	Chrysene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
37	Cyanide	Extraction, Distillation, Colorimetric Method ^[27,28,29]
38	2,4-D	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
39	DDD	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
40	DDE	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
41	DDT	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
42	Dibenz(a,h)anthracene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
43	Di-n-Butyl Phthalate	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
44	1,2-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
45	1,3-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
46	1,4-Dichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
47	3,3-Dichlorobenzidine	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
48	1,1-Dichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
49	1,2-Dichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
50	1,1-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
51	cis-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
52	trans-1,2-Dichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
53	2,4-Dichlorophenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
54	1,2-Dichloropropane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
55	1,3-Dichloropropane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
56	1,3-Dichloropropene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
57	Dieldrin	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
58	Diethyl Phthalate	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
59	2,4-Dimethylphenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
60	2,4-Dinitrophenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
61	2,4-Dinitrotoluene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
62	2,6-Dinitrotoluene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
63	Di-n-Octyl Phthalate	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
64	Endosulfan	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
65	Endrin	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
66	Ethylbenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
67	Fluoranthene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
68	Fluorene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
69	Heptachlor	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
70	Heptachlor epoxide	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
71	Hexachlorobenzene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
72	Hexachloro-1,3-butadiene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
73	n-Hexane	1) Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25] 2) Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method ^[13]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
74	α -HCH	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
75	β -HCH	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
76	γ -HCH	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
77	Hexachlorocyclopentadiene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
78	Hexachloroethane	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
79	Indeno(1,2,3-cd)pyrene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
80	Isophorone	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
81	Lead	1) Digestion, Inductively Coupled Plasma Method ^[7,16] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]
82	Manganese	1) Digestion, Inductively Coupled Plasma Method ^[7,16] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]
83	Mercury	1) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[20] 2) Thermal Decomposition, Amalgamation, and Atomic Absorption Spectrophotometry ^[21] 3) Digestion, Cold-Vapor Atomic Fluorescence Spectrometric Method ^[30]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
84	Methanol	1) Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25] 2) Equilibrium Headspace, Gas Chromatographic/ Mass Spectrometric Method ^[13,25]
85	Methoxychlor	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
86	Methyl Bromide	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
87	Methylene Chloride	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
88	2-methylphenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
89	2-Methylnaphthalene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
90	Methyl tert-Butyl Ether	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
91	Naphthalene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
92	Nickel	1) Digestion, Inductively Coupled Plasma Method ^[7,16] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]
93	Nitrobenzene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
94	N-Nitrosodiphenylamine	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
95	N-Nitrosodi-n-propylamine	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
96	Polychlorinated biphenyls (PCBs) - Aroclor 1016 - Aroclor 1221 - Aroclor 1232 - Aroclor 1242 - Aroclor 1248 - Aroclor 1254 - Aroclor 1260 - 2-Chlorobiphenyl - 2,2',3,5'-Tetrachlorobiphenyl - 2,2',5,5'-Tetrachlorobiphenyl - 2,3',4,4'-Tetrachlorobiphenyl - 2,2',3,4,5'-Pentachlorobiphenyl - 2,2',4,5,5'-Pentachlorobiphenyl - 2,3,3',4',6-Pentachlorobiphenyl - 2,2',3,4,4',5'-Hexachlorobiphenyl - 2,2',3,4,5,5'-Hexachlorobiphenyl - 2,2',3,5,5',6-Hexachlorobiphenyl - 2,2',4,4',5,5'-Hexachlorobiphenyl - 2,2',3,3',4,4',5-Heptachlorobiphenyl - 2,2',3,4,4',5,5'-Heptachlorobiphenyl - 2,2',3,4,4',5',6-Heptachlorobiphenyl - 2,2',3,4',5,5',6-Heptachlorobiphenyl - 2,2',3,3',4,4',5,5',6-Nonachlorobiphenyl	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
97	Pentachlorophenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
98	Phenanthrene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
99	Phenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
100	Pyrene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
101	Selenium	1) Digestion, Inductively Coupled Plasma Method ^[7,16] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]
102	Silver	1) Digestion, Inductively Coupled Plasma Method ^[7,16] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]
103	Styrene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
104	1,1,2,2-Tetrachloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
105	Tetrachloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
106	Toluene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
107	Toxaphene	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
108	TPH (C ₅ -C ₈)	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
109	TPH (C ₈ - C ₁₆)	1) Automate Extraction, Gas Chromatographic Method ^[11,22] 2) Solvent Extraction, Gas Chromatographic Method ^[12,22] 3) Ultrasonic Extraction, Gas Chromatographic Method ^[22,31]
110	TPH (C ₁₆ - C ₃₅)	1) Automate Extraction, Gas Chromatographic Method ^[11,22] 2) Solvent Extraction, Gas Chromatographic Method ^[12,22] 3) Ultrasonic Extraction, Gas Chromatographic Method ^[22,31]
111	1,2,4-Trichlorobenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
112	1,1,1-Trichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
113	1,1,2-Trichloroethane	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
114	Trichloroethylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
115	2,4,5-Trichlorophenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
116	2,4,6-Trichlorophenol	1) Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[10,26] 2) Automated Soxhlet Extraction, Gas Chromatographic/ Mass Spectrometric Method ^[11,26]
117	1,3,5-Trimethylbenzene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
118	Vanadium	1) Digestion, Inductively Coupled Plasma Method ^[7,16] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]
119	Vinyl Acetate	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
120	Vinyl Chloride	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
121	m-Xylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
122	o-Xylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
123	p-Xylene	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
124	Xylene (Total)	Purge and Trap, Gas Chromatographic/ Mass Spectrometric Method ^[15,25]
125	Zinc	1) Digestion, Inductively Coupled Plasma Method ^[7,16] 2) Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[7,17]

เอกสารอ้างอิง

- กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2566. เรื่อง การจัดการสิ่งปฏิกูลหรือวัสดุที่ไม่ใช้แล้ว. ราชกิจจานุเบกษา. 31 พฤษภาคม 2566. เล่มที่ 140 ตอนพิเศษ 126 ง.
- กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2549. เรื่อง กำหนดค่าปริมาณเขม่าควันที่เจือปนในอากาศที่ระบายออกจากปล่องของหม้อน้ำโรงสีข้าวที่ใช้กลบเป็นเชื้อเพลิง. ราชกิจจานุเบกษา. 4 ธันวาคม 2549. เล่มที่ 123 ตอนพิเศษ 125 ง.
- สมาคมวิศวกรรมสิ่งแวดล้อมแห่งประเทศไทย. คู่มือวิเคราะห์น้ำเสีย. พิมพ์ครั้งที่ 4. กรุงเทพฯ: เรือนแก้วการพิมพ์, 2547.
- APHA, AWWA, WEF. Standard Methods for the Examination of Water and Wastewater. 24th ed. Washington, DC: APHA, 2023.

5. United States Environmental Protection Agency. **Standards of Performance for New Stationary Sources**. 40 CFR 60. Appendix A, 2023.

6. United States Environmental Protection Agency. **Test Methods for Evaluating Solid Waste Physical/Chemical Methods**. SW-846, 2014.

7. United States Environmental Protection Agency. **Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Acid Digestion of Sludges and Sediments and Soils**. SW-846 Method 3050B, 1996.

8. United States Environmental Protection Agency. **Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Alkaline Digestion for Hexavalent Chromium**. SW-846 Method 3060A, 1996.

9. United States Environmental Protection Agency. **Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Separatory Funnel Liquid-Liquid Extraction**. SW-846 Method 3510C, 1996.

10. United States Environmental Protection Agency. **Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Soxhlet Extraction**. SW-846 Method 3540C, 1996.

11. United States Environmental Protection Agency. **Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Automated Soxhlet Extraction**. SW-846 Method 3541, 1994.

12. United States Environmental Protection Agency. **Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Microscale Solvent Extraction (MSE)**. SW-846 Method 3570, 2002.

13. United States Environmental Protection Agency. **Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Volatile Organic Compounds (VOCs) in Various Sample Matrices Using Equilibrium Headspace Analysis**. SW-846 Method 5021A, 2014.

14. United States Environmental Protection Agency. **Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Purge-and-Trap for Aqueous Samples**. SW-846 Method 5030B, 1996.

15. United States Environmental Protection Agency. **Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Closed-System Purge-and-Trap and Extraction for Volatile Organics in Soil and Waste Samples**. SW-846 Method 5035, 1996.

16. United States Environmental Protection Agency. **Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Inductively Coupled Plasma-Atomic Emission Spectrometry**. SW-846 Method 6010B, 1996.

17. United States Environmental Protection Agency. **Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Inductively Coupled Plasma-Mass Spectrometry**. SW-846 Method 6020A, 2007.

18. United States Environmental Protection Agency. **Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Antimony and Arsenic (Atomic Absorption, Borohydride Reduction)**. SW-846 Method 7062, 1994. เพิ่มใหม่

19. United States Environmental Protection Agency. **Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Chromium, Hexavalent (Colorimetric)**. SW-846 Method 7196A, 1992.

20. United States...

20. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Mercury in Solid or Semisolid Waste (Manual Cold-Vapor Technique). SW-846 Method 7471B, 2007.

21. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Mercury in Solids and Solutions by Thermal Decomposition, Amalgamation, and Atomic Absorption Spectrophotometry. SW-846 Method 7473, 2007.

22. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Nonhalogenated Organics by Gas Chromatography. SW-846 Method 8015C, 2007.

23. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. pH Electrometric Measurement. SW-846 Method 9040C, 2004.

24. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Soil and Waste pH. SW-846 Method 9045D, 2004.

25. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Volatile Organic Compounds by Gas Chromatography/Mass Spectrometry (GC/MS). SW-846 Method 8260D, 2018.

26. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Semivolatile Organic Compounds by Gas Chromatography/Mass Spectrometry (GC/MS). SW-846 Method 8270E, 2018.

27. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Total and Amenable Cyanide: Distillation SW-846 Method 9010B, 1996.

28. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Cyanide Extraction Procedure for Solids and Oil. SW-846 Method 9013A, 1996.

29. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Cyanide in Waters and Extracts Using Titrimetric and Manual Spectrophotometric Procedures. SW-846 Method 9014, 2014.

30. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Mercury in Sediment and Tissue Samples by Atomic Fluorescence Spectrometry. SW-846 Method 7474, 2007.

31. United States Environmental Protection Agency. Test Methods for Evaluating Solid Waste Physical/Chemical Methods. Ultrasonic Extraction. SW-846 Method 3550C, 2007.





ที่ อก ๐๓๑๐(๑)/ ๔๑๒๑

กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๒๕ เมษายน ๒๕๖๓

เรื่อง เปลี่ยนแปลงบุคลากรของห้องปฏิบัติการวิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๒๙ มีนาคม ๒๕๖๓

ตามคำขอที่อ้างถึง บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด ห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ว-๒๐๔ สถานที่ตั้งเลขที่ ๑๐๔ ซอยพัฒนาการ ๔๐ ถนนพัฒนาการ แขวงพัฒนาการ เขตสวนหลวง กรุงเทพมหานคร ขอเปลี่ยนแปลงบุคลากร ความละเอียดแจ้งแล้ว นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว มีความเห็นดังนี้

๑. ให้ยกเลิกเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๓ ราย

- | | |
|--------------------------|----------------------------|
| ๑) นางสาวพรรณธิดา พุ่มคง | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๕ |
| ๒) นายกำชัย สุทธะ | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๒๑ |
| ๓) นางสาวศุภรดา บัณมยุรา | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๓๘ |

๒. ให้เพิ่มเจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกชน จำนวน ๑๒ ราย

- | | |
|-----------------------------|----------------------------|
| ๑) นางสาวฐานิดา กลิ่นเขียว | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๒ |
| ๒) นางสาวกัญญภััสสร สายคำ | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๓ |
| ๓) นางสาวณัฐนันท์ กันทะวงศ์ | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๔ |
| ๔) นายอำนาจ วงษาเคน | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๕ |
| ๕) นายกฤษณพล ปัญญาวงศ์ | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๖ |
| ๖) นายณชากร หารธา | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๗ |
| ๗) นายวัชรินทร์ ผ่องสามสวน | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๘ |
| ๘) นายณัฐพงศ์ โสภา | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๙ |
| ๙) นายศักรินทร์ ปานเพ็ง | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๙๐ |
| ๑๐) นายณัฐพล ชุ่มชื่น | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๙๑ |
| ๑๑) นายธนา สุพาพันธุ์ | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๙๒ |
| ๑๒) นายนราธร แก้วพงษ์ชา | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๙๓ |

อนึ่ง หนังสือฉบับนี้จะหมดอายุพร้อมหนังสือต่ออายุรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
ในวันที่ ๒ กันยายน ๒๕๖๙

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ



(นายพริยศ กลั่นกรอง)

รองอธิบดี ปฏิบัติราชการแทน
อธิบดีกรมโรงงานอุตสาหกรรม

กองวิจัยและเตือนภัยมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕

โทรสาร ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๙๙

ไปรษณีย์อิเล็กทรอนิกส์ saraban@diw.mail.go.th



ที่อก ๐๓๑๐(๑)/ ๑๒ ๓๖ ๘ /



กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๑๘ ธันวาคม ๒๕๖๗

เรื่อง ยกเลิกบุคลากรของห้องปฏิบัติการวิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๒ ธันวาคม ๒๕๖๗


ตามคำขอที่อ้างถึง บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด ห้องปฏิบัติการ
วิเคราะห์เอกชน เลขทะเบียน ว-๒๐๔ สถานที่ตั้งเลขที่ ๑๐๔ ซอยพัฒนาการ ๔๐ ถนนพัฒนาการ แขวงพัฒนาการ
เขตสวนหลวง กรุงเทพมหานคร ขอยกเลิกบุคลากร ความละเอียดแจ้งแล้ว นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้ยกเลิกเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์
จำนวน ๘ ราย ได้แก่

- | | |
|-------------------------------|----------------------------|
| ๑) นายประพจน์ วรรณชูชัย | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๖๐ |
| ๒) นายจิรณัฐ ขาวละออ | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๐๗๒ |
| ๓) นายพีรพัฒน์ กำคำ | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๐๘ |
| ๔) นางสาวอรยา คำคล่อง | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๓๔ |
| ๕) นายกิตติพงศ์ แซ่ลี | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๔๔ |
| ๖) นายจิรเมธ ประเสริฐศิริพงศ์ | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๐ |
| ๗) นายภัทรพงษ์ มณฑาทอง | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๖๗ |
| ๘) นางสาวจารุวรรณ กระจำงพันธ์ | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๘๑ |

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

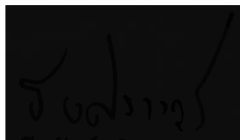

(นายธรรตศน์ อัครางกูร ณ อยุธยา)

รองอธิบดี ปฏิบัติราชการแทน

อธิบดีกรมโรงงานอุตสาหกรรม

กองวิจัยและเตือนภัยมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕

โทรสาร ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๙๙

ไปรษณีย์อิเล็กทรอนิกส์ saraban@diw.mail.go.th



“อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว”



ที่ อก ๐๓๑๐(๑)/ ๓ ๑ ๕ ๐



กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๕๐๐

๑๐ เมษายน ๒๕๖๔

เรื่อง ยกเลิกบุคลากรของห้องปฏิบัติการวิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๒ เมษายน ๒๕๖๔

ตามคำขอที่อ้างถึง บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด ห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ว-๒๐๔ สถานที่ตั้งเลขที่ ๑๐๔ ซอยพัฒนาการ ๔๐ ถนนพัฒนาการ แขวงพัฒนาการ เขตสวนหลวง กรุงเทพมหานคร ขอยกเลิกบุคลากร ความละเอียดแจ้งแล้ว นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้ยกเลิกเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์
จำนวน ๒ ราย ได้แก่

๑) นายธิตินงค์ บัวแดง

ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๐๒

๒) นายมงคล ผลาทิพย์

ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๑๐

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ



(นายธีรทัศน์ อิศรางกูร ณ อยุธยา)

รองอธิบดี ปฏิบัติราชการแทน

อธิบดีกรมโรงงานอุตสาหกรรม

กองวิจัยและเตือนภัยมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕

โทรสาร ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๙๙

ไปรษณีย์อิเล็กทรอนิกส์ saraban@diw.mail.go.th



ที่ อก ๐๓๑๐(๑)/ ๕๖ ๗ ๓



กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๐๕ มิถุนายน ๒๕๖๔

เรื่อง เปลี่ยนแปลงบุคลากรของห้องปฏิบัติการวิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๗ พฤษภาคม ๒๕๖๔

ตามคำขอที่อ้างถึง บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด ห้องปฏิบัติการ
วิเคราะห์เอกชน เลขทะเบียน ว-๒๐๔ สถานที่ตั้งเลขที่ ๑๐๔ ซอยพัฒนาการ ๔๐ ถนนพัฒนาการ แขวงพัฒนาการ
เขตสวนหลวง กรุงเทพมหานคร ขอเปลี่ยนแปลงบุคลากร ความละเอียดแจ้งแล้ว นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้เพิ่มเจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกชน
จำนวน ๑๑ ราย ได้แก่

- | | |
|-----------------------------------|----------------------------|
| ๑) นายคุณากร มั่นชื่น | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๙๔ |
| ๒) นายชัยมงคล แสนมาตร | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๙๕ |
| ๓) นายเอกรินทร์ บุตส์ักดิ์ | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๙๖ |
| ๔) นายพิชานนท์ อินปริง | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๙๗ |
| ๕) นายศตวรรษ แก้วกันหา | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๙๘ |
| ๖) นายวิกรม มีศิริ | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๑๙๙ |
| ๗) นายภคนันท์ คำจันทราช | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๒๐๐ |
| ๘) นายศิริวิทย์ มีไพฑูรย์ | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๒๐๑ |
| ๙) นายธีรพงษ์ ศรีคำแหง | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๒๐๒ |
| ๑๐) นายอภิสิทธิ์ ศรีคงแก้ว | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๒๐๓ |
| ๑๑) ว่าที่ร้อยตรี ภาณุพงศ์ แสนศรี | ทะเบียนเลขที่ ว-๒๐๔-จ-๐๒๐๔ |

อนึ่ง หนังสือฉบับนี้จะสิ้นอายุพร้อมหนังสือต่ออายุรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
ในวันที่ ๒ กันยายน ๒๕๖๔

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ



(นายศิระ จันทรเลิศ)

นักวิทยาศาสตร์เชี่ยวชาญ รัชการราชการแทน
ผู้อำนวยการกองวิจัยและเตือนภัยมลพิษโรงงาน
ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม

กองวิจัยและเตือนภัยมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕

โทรสาร ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๙๙

ไปรษณีย์อิเล็กทรอนิกส์ saraban@diw.mail.go.th



“อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว”





๐๔ สิงหาคม ๒๕๖๗

เรื่อง ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

เรียน กรรมการผู้จัดการ บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์
เอกชน ลงวันที่ ๒๗ พฤษภาคม ๒๕๖๗

สิ่งที่ส่งมาด้วย เอกสารแนบท้ายหนังสือต่ออายุรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด จำนวน ๓ แผ่น

ตามคำขอที่อ้างถึง บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด ขอต่ออายุ
หนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ว-๓๒๓ สถานที่ตั้งเลขที่ ๖๑๖/๑๐ หมู่ที่ ๕
ตำบลแม่น้ำคู้ อำเภอลวกแดง จังหวัดระยอง ต่อกรมโรงงานอุตสาหกรรม นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย)
จำกัด ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน โดยมีองค์ประกอบดังนี้

ก. ผู้ควบคุมห้องปฏิบัติการวิเคราะห์เอกชน

- | | |
|--------------------------|----------------------------|
| ๑) นายเดช ช้างชน | ทะเบียนเลขที่ ว-๓๒๓-ค-๐๐๐๑ |
| ๒) นางวิลาวัลย์ บริรักษ์ | ทะเบียนเลขที่ ว-๓๒๓-ค-๐๐๐๒ |
| ๓) นายสุพจน์ สลามเต๊ะ | ทะเบียนเลขที่ ว-๓๒๓-ค-๐๐๐๓ |

ข. เจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกชน

- | | |
|--------------------------------|----------------------------|
| ๑) นายณัฐพงษ์ เพ็งขาวนา | ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๐๑ |
| ๒) นางสาวกัลยพรรณ รักษ์ดี | ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๐๒ |
| ๓) นางสาวจุฑารัตน์ สิททองกลาง | ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๐๓ |
| ๔) นางสาวจิตสุภา ประเทืองสุข | ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๐๔ |
| ๕) นายสรสรเสรีญ์ ค่อยยกสุย | ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๐๕ |
| ๖) นายณัฐวุฒิ ออมพรมราช | ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๐๖ |
| ๗) นายจิตรกร สีวะสา | ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๐๗ |
| ๘) นายสิทธิปวิญญ์ สุวรรณรัตน์ | ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๐๘ |
| ๙) นายสิทธิพันธ์ เสนาชีว | ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๐๙ |
| ๑๐) นายอนุเวศน์ เตมา | ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๑๐ |
| ๑๑) นายสุรวิทย์ นราพงษ์ | ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๑๑ |
| ๑๒) นายณัฐพล เจียงวรวงศ์ | ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๑๓ |
| ๑๓) นายชานนท์ บุญชื่น | ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๑๔ |
| ๑๔) นายณัฐกานต์ วงศ์อินทร์อยู่ | ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๑๕ |
| ๑๕) นายอานนท์ โพธิ์พระทอง | ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๑๖ |

๑๖) นายณัชพล...

- ๑๖) นายณัฏฐพล ถ้ำกลาง
- ๑๗) นายศุภณัฐ พิสัยพันธ์
- ๑๘) นายวสันต์ คินันติ
- ๑๙) นายวรัญญู นิมพาสิทธิ์
- ๒๐) นายศุภณัฐ สกฤกิตติมงคล
- ๒๑) นายเอกชัย ถิ่นทอง
- ๒๒) นายพงษ์เทพ สิทธิเลาะ
- ๒๓) นายทินกร กุมภาชี
- ๒๔) นางสาวนันทยา เบญจพันธ์
- ๒๕) นายสิทธิชัย ยันพิมาย
- ๒๖) นางสาวปภาณิน หลอดทอง
- ๒๗) นางสาวพจนา สีดา
- ๒๘) นางสาวธนิดา กุลศิริวงศ์
- ๒๙) นายพิทยา ทองแดง
- ๓๐) นางสาวชลธิชา สุปงกช
- ๓๑) ว่าที่ร้อยตรี ธรณชัย ม่วงมา
- ๓๒) นายวรารุณ พับพา
- ๓๓) นายศักดิ์นรินทร์ จรัสกาย
- ๓๔) นายสุรศักดิ์ สาชิน
- ๓๕) นายสถาพร ถาแก้ว
- ๓๖) นายสุทธิดำรง โชคปิตินันท์
- ๓๗) นายวัลลภ หันไชยเนาว์
- ๓๘) นางสาววนาลี เจริญตระกูล
- ๓๙) นายธนะสิทธิ์ วงศ์ไชย
- ๔๐) นายชัยนุสรณ์ เลิศนันทกุลชัย
- ๔๑) นายสัจจา เพ็ชรแสง
- ๔๒) นายกัณตภณ มณีสัมพันธ์
- ๔๓) นายธารินทร์ อ็อกจินดา
- ๔๔) นายศุภชัย วงศ์สุริยฉาย
- ๔๕) นายไสว ตันโพธิ์
- ๔๖) นางสาวกิตติยา สันญาอริยาภรณ์
- ๔๖) นางสาวธิดารัตน์ ศิริมังคะโร
- ๔๗) นายพิพัฒน์ นิภัทร์เศรษฐ์
- ๔๘) นายศิริวิทย์ เรืองสม
- ๔๙) นายปารเมศ สัตยาคุณ
- ๕๐) นายนฤนาท ธรรมสระโร
- ๕๑) นางสาวศุภรัตน์ โสจันทร์

- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๑๗
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๑๘
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๑๙
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๒๐
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๒๑
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๒๒
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๒๓
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๒๔
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๒๕
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๒๖
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๒๗
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๒๘
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๒๙
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๓๐
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๓๑
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๓๒
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๓๓
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๓๔
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๓๕
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๓๖
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๓๗
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๓๘
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๓๙
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๔๐
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๔๑
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๔๒
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๔๓
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๔๔
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๔๕
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๔๖
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๔๗
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๔๘
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๔๙
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๕๐
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๕๑
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๕๒
- ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๕๓

๕๒) นายพชรกร เจ็จเจริญ	ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๕๔
๕๓) นายทิวากร เชื้อมาก	ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๕๕
๕๔) นายอนุรักษ์ ทองขจรศักดิ์	ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๕๖
๕๕) นายอภิชาติ วิชาส	ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๕๗
๕๖) นายจรัสระวี ศรีรักษา	ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๕๘
๕๘) นายประสานมิตร เชื้อนเพชร	ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๕๙
๕๙) นายภาณุวัฒน์ วังบง	ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๖๐
๖๐) นายสันติ ชัยชนะ	ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๖๑
๖๑) นายทินกร กุลชาติ	ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๖๒

ค. ขอบข่ายชนิดสารมลพิษที่ได้รับขึ้นทะเบียนให้วิเคราะห์ในน้ำเสีย น้ำใต้ดิน อากาศเสีย ตามสิ่งที่ส่งมาด้วย

หนังสือฉบับนี้จะหมดอายุในวันที่ ๒๘ มิถุนายน ๒๕๗๑ หากประสงค์จะต่ออายุหนังสือ รับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ให้ยื่นคำขอต่ออายุพร้อมเอกสารประกอบคำขอต่อกรมโรงงาน อุตสาหกรรมภายใน ๖๐ วัน ก่อนวันสิ้นอายุของหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ



(นายพริต กุลกรอง)

รองอธิบดี ปฏิบัติราชการแทน

อธิบดีกรมโรงงานอุตสาหกรรม

ศูนย์วิจัยและเตือนภัยมลพิษโรงงานภาคตะวันออก

โทร. ๐ ๓๓๑๓ ๖๐๕๔ ต่อ ๕๐๐๑-๒

ไปรษณีย์อิเล็กทรอนิกส์ eirw@diw.mail.go.th



“อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว”



เอกสารแนบท้ายหนังสือเปลี่ยนแปลงสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด เลขทะเบียน ว-๓๒๓
ที่ ออก ๐๓๒๐/ ๗ ๙๓ ๘ ลงวันที่ ๐๘ สิงหาคม ๒๕๖๗

ขอขยายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๒๔ รายการ
น้ำเสีย จำนวน 14 รายการ

ลำดับ ที่	สารมลพิษ	วิธีวิเคราะห์
1	Biochemical Oxygen Demand	1) 5-Day BOD Test, Membrane Electrode Method ^[2] 2) 5-Day BOD Test, Azide Modification Method ^[2]
2	Chemical Oxygen Demand	1) Open Reflux, Titrimetric Method ^[2] 2) Closed Reflux, Colorimetric Method ^[2] 3) Closed Reflux, Titrimetric Method ^[2]
3	Color	ADMI Weighted-Ordinate Spectrophotometric Method ^[2]
4	Cyanide	Distillation, Colorimetric Method ^[2]
5	Formaldehyde	Distillation, Colorimetric Method ^[1]
6	Free Chlorine	DPD Ferrous Titrimetric Method ^[2]
7	Oil and Grease	Liquid-Liquid, Partition-Gravimetric Method ^[2]
8	pH	Electrometric Method ^[2]
9	Phenols	1) Distillation, Chloroform Extraction Method ^[2] 2) Distillation, Direct Photometric Method ^[2]
10	Sulfide	ZnS Precipitation, Iodometric Method ^[2]
11	Temperature	Field Method ^[2]
12	Total Dissolved Solids	Dried at 180 °C ^[2]
13	Total Kjeldahl Nitrogen	Semi-Macro Kjeldahl Method ^[2]
14	Total Suspended Solids	Dried at 103-105 °C ^[2]

น้ำใต้ดิน จำนวน 3 รายการ

ลำดับ ที่	สารมลพิษ	วิธีวิเคราะห์
1	Cyanide	Distillation, Colorimetric Method ^[2]
2	pH	Electrometric Method ^[2]
3	Phenols	Distillation, Direct Photometric Method ^[2]

อากาศเสีย...

อากาศเสีย (ปล่อยระบาย) จำนวน 7 รายการ

ลำดับ ที่	สารมลพิษ	วิธีวิเคราะห์
1	Carbon Monoxide	1) Sampling Bag, Non-Dispersive Infrared Method ^[5] 2) Instrumental Analyzer Method ^[9]
2	Hydrogen Sulfide	Absorption Sampling, Iodometric Method ^[5]
3	Opacity	Ringelmann's Method ^[3,4]
4	Oxide of Nitrogen	1) Absorption Sampling, Phenoldisulfonic Acid Method ^[8] 2) Instrumental Analyzer Method ^[10]
5	Sulfur Dioxide	1) Absorption Sampling, Barium-Thorin Titrimetric Acid Method ^[5] 2) Instrumental Analyzer Method ^[11]
6	Sulfuric Acid	Isokinetic Sampling, Barium – Titrimetric Method ^[6]
7	Total Suspended Particulate	Isokinetic Sampling, Gravimetric Method ^[7]

เอกสารอ้างอิง

1. รัชชัย พรหมสวัสดิ์ และวิบูลย์ลักษณ์ วิสุมิศักดิ์, บรรณาธิการ. (2547) คู่มือวิเคราะห์น้ำเสีย. พิมพ์ครั้งที่ 4. กรุงเทพฯ: สมาคมวิศวกรรมสิ่งแวดล้อมแห่งประเทศไทย.

2. APHA, AWWA, WEF. Standard Methods for the Examination of Water and Wastewater. 24th ed. Washington, DC : APHA, 2023

3. กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2549. เรื่อง กำหนดค่าปริมาณเขม่าควันที่เจือปนในอากาศที่ระบายออกจากปล่องของหม้อน้ำโรงสีข้าวที่ใช้กลบเป็นเชื้อเพลิง. ราชกิจจานุเบกษา. 4 ธันวาคม 2549. เล่มที่ 123 ตอนพิเศษ 125ง

4. กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2549. เรื่อง กำหนดค่าปริมาณเขม่าควันที่เจือปนในอากาศที่ระบายออกจากปล่องของหม้อน้ำของโรงงาน. ราชกิจจานุเบกษา. 4 ธันวาคม 2549. เล่มที่ 123 ตอนพิเศษ 125ง

5. United States Environmental Protection Agency. Standards of Performance for New Stationary Sources. 40 CFR 60. Appendix A, 2017.

6. United States Environmental Protection Agency. Standards of Performance for New Stationary Sources. 40 CFR 60. Appendix A, 2019.

7. United States...

7. United States Environmental Protection Agency. **Standards of Performance for New Stationary Sources**. 40 CFR 60. Appendix A, 2020.

8. United States Environmental Protection Agency. **Standards of Performance for New Stationary Sources**. 40 CFR 60. Appendix A, 2023.

9. United States Environmental Protection Agency. **Determination of Carbon Monoxide Emission from Stationary Sources; Instrumental Analyzer Procedure**. 40 CFR 60. Appendix A Method 10, 2017.

10. United States Environmental Protection Agency. **Determination of Oxide of Nitrogen Emission from Stationary Sources; Instrumental Analyzer Procedure**. 40 CFR 60. Appendix A Method 7E, 2023.

11. United States Environmental Protection Agency. **Determination of Sulfur dioxide Emission from Stationary Sources; Instrumental Analyzer Procedure**. 40 CFR 60. Appendix A Method 6C, 2017.

A handwritten signature in black ink, appearing to be 'D. E.', is located in the bottom right corner of the page.

ที่ อก ๐๓๒๐/ ๑๐๐๙ ๙



กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๐๔ ตุลาคม ๒๕๖๗

เรื่อง แกไขรายชื่อเจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกชน

เรียน กรรมการผู้จัดการ บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด

อ้างถึง หนังสือ บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด เลขที่ Env 2024/005

ลงวันที่ ๓๐ สิงหาคม ๒๕๖๗

ตามหนังสือที่อ้างถึง บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด ห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ว-๓๒๓ สถานที่ตั้งเลขที่ ๖๑๖/๑๐ หมู่ที่ ๕ ตำบลแม่น้ำคู้ อำเภอบลวกแดง จังหวัดระยอง ขอแก้ไขชื่อเจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกชน เนื่องจากมีความคลาดเคลื่อน ความละเอียดแจ้งแล้ว นั้น

กรมโรงงานอุตสาหกรรม ได้รับทราบและดำเนินการแก้ไขรายชื่อเจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกชน จำนวน ๕ ราย ตามที่แจ้งเรียบร้อยแล้ว เป็นดังนี้

ลำดับที่ ๒๗ นางพจนา สีดา

ลำดับที่ ๒๘ นางสาวธนิศา กุลสุริวงศ์

ลำดับที่ ๓๐ นางชลธิชา สุปงกช

ลำดับที่ ๓๖ นายสุทธิดำรงค์ โชคปิตินันท์

ลำดับที่ ๔๒ นายกันตภณ มณีสัมพันธ์

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ



(นายพรยศ กลืนกรอง)

รองอธิบดี ปฏิบัติราชการแทน

อธิบดีกรมโรงงานอุตสาหกรรม

ศูนย์วิจัยและเตือนภัยมลพิษโรงงานภาคตะวันออก

โทร. ๐ ๓๓๑๓ ๖๐๕๙ ต่อ ๕๐๐๑-๒

ไปรษณีย์อิเล็กทรอนิกส์ eirw@diw.mail.go.th



“อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว”



ที่ อก ๐๓๑๐(๓)/ ๕๒๔๖



กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๒๐ พฤษภาคม ๒๕๖๔

เรื่อง ยกเลิกบุคลากรของห้องปฏิบัติการวิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๑๐ เมษายน ๒๕๖๔

ตามคำขอที่อ้างถึง บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด ห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ว-๓๒๓ สถานที่ตั้งเลขที่ ๖๑๖/๑๐ หมู่ที่ ๕ ตำบลแม่น้ำคู้ อำเภอลำลูกเกด จังหวัดระยอง ขอยกเลิกบุคลากร ความละเอียดแจ้งแล้ว นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้ยกเลิกเจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกชน จำนวน ๑ ราย ได้แก่ นายปารามาศ สัตยาคุณ ทะเบียนเลขที่ ว-๓๒๓-จ-๐๐๕๑

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

(นายประสม ดำรงพงษ์)

ผู้อำนวยการกองวิจัยและเตือนภัยมลพิษโรงงาน
ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม

ศูนย์วิจัยและเตือนภัยมลพิษโรงงานภาคตะวันออก

โทร. ๐ ๓๓๑๓ ๖๐๕๔ ต่อ ๕๐๐๑-๒

ไปรษณีย์อิเล็กทรอนิกส์ eirw@diw.mail.go.th



“อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว”



ที่ อก ๐๓๑๐(๓)/ ๕๔๐๙



กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๒๗ พฤษภาคม ๒๕๖๘

เรื่อง เปลี่ยนแปลงชื่อ-สกุลบุคลากรของห้องปฏิบัติการวิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๑๕ พฤษภาคม ๒๕๖๘

ตามคำขอที่อ้างถึง บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด ห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ว-๓๒๓ สถานที่ตั้งเลขที่ ๖๑๖/๑๐ หมู่ที่ ๕ ตำบลแม่წყ่า อำเภอปลวกแดง จังหวัดระยอง ขอเปลี่ยนแปลงชื่อ-สกุลบุคลากร ความละเอียดแจ้งแล้ว นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้เปลี่ยนแปลงชื่อ-สกุลบุคลากร จำนวน ๑ ราย
จากนายธนสิทธิ์ วงศ์ไชย เป็น นายอมลวิษฐ์ วงศ์ไชย

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

(นายประสม ดำรงพงษ์)

ผู้อำนวยการกองวิจัยและเตือนภัยมลพิษโรงงาน
ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม

ศูนย์วิจัยและเตือนภัยมลพิษโรงงานภาคตะวันออก

โทร. ๐ ๓๓๑๓ ๖๐๕๙ ต่อ ๕๐๐๑-๒

ไปรษณีย์อิเล็กทรอนิกส์ eirw@diw.mail.go.th



“อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว”



ที่ อก ๐๓๒๒/๑๓๖๕๙



กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๒๕ ก.ย. ๒๕๖๖

เรื่อง ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

เรียน กรรมการผู้จัดการ บริษัท เอแอลเอส แลบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๑๔ กรกฎาคม ๒๕๖๖

สิ่งที่ส่งมาด้วย เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
บริษัท เอแอลเอส แลบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด จำนวน ๓ แผ่น

ตามหนังสือที่อ้างถึง บริษัท เอแอลเอส แลบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด ขอต่ออายุ
หนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ว-๒๖๗ สถานที่ตั้ง เลขที่ ๑๑๔/๑ หมู่ที่ ๘
ถนนกาญจนวนิช ตำบลบ้านพรุ อำเภอหาดใหญ่ จังหวัดสงขลา ต่อกรมโรงงานอุตสาหกรรม นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้บริษัท เอแอลเอส แลบอราทอรี กรุ๊ป (ประเทศไทย)
จำกัด ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน โดยมีองค์ประกอบดังนี้

ก. ผู้ควบคุมดูแลห้องปฏิบัติการวิเคราะห์

นางสาวกนิษฐา เหมประสาทร

ทะเบียนเลขที่ ว-๒๖๗-ค-๐๐๐๑

ข. เจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์

๑) นางสาวอินทิรา คงประยูร

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๐๑

๒) นางสาวอมรรัตน์ เพชรประดับ

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๐๒

๓) นายทักษิณ อินโตรม

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๐๓

๔) นางสาวอณัฏา บุญเพชร

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๐๔

๕) นางสาวสุทธริกา ทิพย์รัตน์

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๐๕

๖) นางสาวนริสา นฤมิตร

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๐๖

๗) นายวุฒิชัย ทวยเจริญ

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๐๗

๘) นายยงศิลป์ รังษี

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๐๘

๙) นายอภิวัฒน์ ฉันทะ

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๐๙

๑๐) นายศิริชัย เกลี้ยงเกิด

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๑๐

๑๑) นายสมศักดิ์ จันทรวง

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๑๑

๑๒) นางสาวพิชญา ศุภรานนท์

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๑๒

๑๓) นายปัญญา เกียรติพิมุรักษ์

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๑๓

๑๔) นางสาวศศิณีภา รอดทองอ่อน

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๑๔

๑๕) นางสาวชุติมา สุขสวัสดิ์

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๑๕

๑๖) นางสาวจันทิมา คงทน

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๑๖

๑๗) นางสาวกุลวดี เรืองประพันธ์

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๑๗

๑๘) นางสาวอาทิตย์ยา เมืองแก้ว

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๑๘

๑๙) นางสาวกวิณณา ฉุนย่อง

ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๑๙

ค. ขอบข่ายสารมลพิษที่ได้รับขึ้นทะเบียนให้วิเคราะห์ในน้ำเสียและอากาศเสีย ตามสิ่งที่ส่งมาด้วย

หนังสือฉบับนี้....



หนังสือฉบับนี้จะหมดอายุในวันที่ ๑๖ สิงหาคม ๒๕๖๙ หากประสงค์จะต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ให้ยื่นคำขอต่ออายุพร้อมเอกสารประกอบคำขอต่อกรมโรงงานอุตสาหกรรม ภายใน ๓๐ วัน ก่อนวันสิ้นอายุของหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ทั้งนี้ สามารถยื่นคำขอผ่านระบบอิเล็กทรอนิกส์ได้ที่หน้าเว็บไซต์กรมโรงงานอุตสาหกรรม

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ



(นายเนเรศวร์ ตริยงค์)

ผู้อำนวยการศูนย์วิจัยและเตือนภัยมลพิษโรงงานภาคใต้
ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม

ศูนย์วิจัยและเตือนภัยมลพิษโรงงานภาคใต้
โทร. ๐ ๙๔๓๒ ๕๐๒๙, ๐ ๙๔๘๙ ๐๖๓๔ ต่อ ๕๒๐๑
ไปรษณีย์อิเล็กทรอนิกส์ sirw@diw.mail.go.th

เอกสารแนบท้ายหนังสือรับต่ออายุขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด เลขทะเบียน ว-๒๖๗
ที่ อก ๐๓๒๒/๑๓๖๕๔ ลงวันที่ ๒๕ ก.ย. ๒๕๖๖

ขอขยายสารมลพิษที่ได้รับการขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๓๗ รายการ
น้ำเสีย จำนวน 25 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Arsenic	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1]
2	Barium	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1]
3	Biochemical Oxygen Demand	5-Day BOD Test, Azide Modification Method ^[1] 5-Day BOD Test, Membrane Electrode Method ^[1]
4	Cadmium	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1]
5	Chemical Oxygen Demand	Closed Reflux, Colorimetric Method ^[1] Closed Reflux, Titrimetric Method ^[1]
6	Chromium	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1]
7	Color	ADMI Weighted-Ordinate Spectrophotometric Method ^[1]
8	Copper	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1]
9	Formaldehyde	Distillation, Colorimetric Method ^[2]
10	Free Chlorine	DPD Ferrous Titrimetric Method ^[1]
11	Hexavalent Chromium	Filtration, Colorimetric Method ^[1]
12	Lead	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1]
13	Manganese	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1]
14	Mercury	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1]
15	Nickel	Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[1]
16	Oil & Grease	Liquid-Liquid, Partition-Gravimetric Method ^[1]

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
17	pH	Electrometric Method ^[1]
18	Phenol	Distillation, Direct Photometric Method ^[1]
19	Selenium	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[1]
20	Sulfide	ZnS Precipitation, Iodometric Method ^[1]
21	Temperature	Laboratory and Field Methods ^[1]
22	Total Dissolved Solids	Dried at 180 °C ^[1]
23	Total Suspended Solids	Dried at 103-105 °C ^[1]
24	Trivalent Chromium	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[1]
25	Zinc	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method ^[1]

อากาศเสีย จำนวน 12 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Antimony	Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[3]
2	Arsenic	Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[3]
3	Carbon Monoxide	Sampling Bag Non-Dispersive Infrared Method ^[3]
4	Copper	Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[3]
5	Dioxins	Isokinetic Sampling, Analysis by ISO/IEC 17025 Accredited Laboratory or Analysis by Department of Industrial Works Registered Laboratory ^[3]
6	Hydrogen Sulfide	Absorption, Iodometric Method ^[3]
7	Lead	Isokinetic Sampling, Digestion, Inductively Coupled Plasma/Mass Spectrometric Method ^[3]
8	Opacity	Ringelmann's Method ^[4]
9	Oxides of Nitrogen	Absorption Sampling, Phenoldisulfonic acid Method ^[3]
10	Sulfur Dioxide	Absorption Sampling, Barium-Thorin Titrimetric Method ^[3]
11	Sulfuric acid	Isokinetic Sampling, Barium-Thorin Titrimetric Method ^[3]
12	Total Suspended Particulate	Isokinetic Sampling, Gravimetric Method ^[3]

เอกสารอ้างอิง

1. APHA, AWWA, WEF. Standard Methods for the Examination of Water and Wastewater. 24th ed. Washington, DC: APHA, 2023.
2. สมาคมวิศวกรรมสิ่งแวดล้อมแห่งประเทศไทย. คู่มือวิเคราะห์น้ำเสีย. พิมพ์ครั้งที่ 4. กรุงเทพฯ: เรือนแก้วการพิมพ์, 2547.
3. United States Environmental Protection Agency. Standards of Performance for New Stationary Sources. 40 CFR 60. Appendix A, 2020.
4. กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2549. เรื่อง กำหนดค่าปริมาณเขม่าควันที่เจือปนในอากาศที่ระบายออกจากปล่องของหม้อน้ำโรงสีข้าวที่ใช้แกลบเป็นเชื้อเพลิง. ราชกิจจานุเบกษา. 4 ธันวาคม 2549. เล่มที่ 123 ตอนพิเศษ 125ง.


(นางสาวบุษยา รัตนสุภา)
นักวิทยาศาสตร์ชำนาญการ



ที่ อก ๐๓๒๒/๔๐๐๙ .1

กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๒ ๓ เมษายน ๒๕๖๗

เรื่อง เปลี่ยนแปลงบุคลากรของห้องปฏิบัติการวิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด

อ้างถึง คำขอเปลี่ยนแปลงบุคลากรของห้องปฏิบัติการวิเคราะห์เอกชน ลงวันที่ ๒๙ มีนาคม ๒๕๖๗

ตามคำขอฯ ที่อ้างถึง บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด ห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ว-๒๖๗ สถานที่ตั้งเลขที่ ๑๑๔/๑ หมู่ที่ ๘ ถนนกาญจนวนิช ตำบลบ้านพรุ อำเภอกาบังใหญ่ จังหวัดสงขลา ขอเปลี่ยนแปลงบุคลากรของห้องปฏิบัติการวิเคราะห์ ความละเอียดแจ้งแล้ว นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว เห็นชอบให้เพิ่มเจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกชน จำนวน ๓ ราย ได้แก่

- | | |
|--------------------------|----------------------------|
| ๑) นายฟูรคอน เกษตรกาลาม์ | ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๒๐ |
| ๒) นายพลเทพ สืบแก้ว | ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๒๑ |
| ๓) นายธนาวุฒิ ปิ่นทอง | ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๒๒ |

อนึ่ง หนังสือเห็นชอบฉบับนี้จะสิ้นอายุพร้อมหนังสือขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชนในวันที่ ๑๖ สิงหาคม ๒๕๖๙ ทั้งนี้ หากประสงค์จะต่ออายุ สามารถยื่นคำขอผ่านระบบอิเล็กทรอนิกส์ได้ที่หน้าเว็บไซต์กรมโรงงานอุตสาหกรรม

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ



(นายพรยศ กลิ่นกรอง)

รองอธิบดี ปฏิบัติราชการแทน

อธิบดีกรมโรงงานอุตสาหกรรม

ศูนย์วิจัยและเตือนภัยมลพิษโรงงานภาคใต้

โทร. ๐ ๗๔๓๒ ๕๐๒๙, ๐ ๗๔๘๙ ๐๖๓๔ ต่อ ๕๒๐๑

ไปรษณีย์อิเล็กทรอนิกส์ sirw@diw.mail.go.th



ที่ ออก ๐๓๑๐(๕)/ ๑ ๓ ๒



กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๐๗ มกราคม ๒๕๖๘

เรื่อง เปลี่ยนแปลงบุคลากรและสารมลพิษที่วิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท เอแอลเอส แลบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๒ ธันวาคม ๒๕๖๗

สิ่งที่ส่งมาด้วย เอกสารแนบท้ายหนังสือเปลี่ยนแปลงบุคลากรและสารมลพิษที่วิเคราะห์
บริษัท เอแอลเอส แลบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด จำนวน ๒ แผ่น

ตามคำขอ ที่อ้างถึง บริษัท เอแอลเอส แลบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด ห้องปฏิบัติการ
วิเคราะห์เอกชน เลขทะเบียน ว-๒๖๗ สถานที่ตั้งเลขที่ ๑๑๔/๑ หมู่ที่ ๘ ถนนกาญจนวนิช ตำบลบ้านพรุ
อำเภอหาดใหญ่ จังหวัดสงขลา ขอเปลี่ยนแปลงบุคลากรและสารมลพิษที่วิเคราะห์ ต่อกรมโรงงานอุตสาหกรรม นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว มีความเห็นดังนี้

๑. ให้ยกเลิกเจ้าหน้าที่ประจำห้องปฏิบัติการวิเคราะห์ จำนวน ๑ ราย ได้แก่

นางสาวอินทรา คงประยูร ทะเบียนเลขที่ ว-๒๖๗-จ-๐๐๐๑

๒. ให้เพิ่มขอบข่ายชนิดสารมลพิษที่วิเคราะห์ในน้ำ/น้ำเสีย และน้ำใต้ดิน ตามสิ่งที่ส่งมาด้วย

อนึ่ง หนังสือฉบับนี้จะสิ้นสุดอายุพร้อมหนังสือต่ออายุรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
ในวันที่ ๑๖ สิงหาคม ๒๕๖๙

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

(นายธีรทัศน์ อิศรางกูร ณ อยุธยา)

รองอธิบดี ปฏิบัติราชการแทน

อธิบดีกรมโรงงานอุตสาหกรรม

กองวิจัยและเตือนภัยมลพิษโรงงาน

ศูนย์วิจัยและเตือนภัยมลพิษโรงงานภาคใต้

โทร. ๐ ๗๔๓๒ ๕๐๒๙, ๐ ๗๔๘๙ ๐๖๓๔ ต่อ ๕๒๐๑

ไปรษณีย์อิเล็กทรอนิกส์ sirw@diw.mail.go.th



เอกสารแนบท้ายหนังสือเปลี่ยนแปลงบุคลากรและสารมลพิษที่วิเคราะห์

บริษัท เอแอลเอส แลบบอราทอรี กรุ๊ป (ประเทศไทย) จำกัด

เลขทะเบียน ว-๒๖๗

ที่ อก ๐๓๑๐(๕)/' ๑ ๓ ๒

ลงวันที่ ๐๗ มกราคม ๒๕๖๘

ขอขยายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๑๙ รายการ

น้ำ/น้ำเสีย จำนวน 1 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Total Kjeldahl Nitrogen	Semi-Micro-Kjeldahl Methods

น้ำใต้ดิน จำนวน 18 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Antimony	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method
2	Arsenic	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method
3	Barium	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method
4	Beryllium	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method
5	Cadmium	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method
6	Chromium	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method
7	Chromium (III)	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method; Colorimetric Method; Calculation
8	Chromium (VI)	Filtration, Colorimetric Method
9	Lead	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method
10	Manganese	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method
11	Mercury	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method
12	Nickel	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
13	pH	Electrometric Method
14	Phenol	Distillation, Direct Photometric Method
15	Selenium	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method
16	Silver	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method
17	Vanadium	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method
18	Zinc	Digestion, Inductively Coupled Plasma/ Mass Spectrometric Method

เอกสารอ้างอิง

APHA, AWWA, WEF. Standard Methods for the Examination of Water and Wastewater. 24th ed. Washington, DC: APHA, 2023.


(นางสาวบุษยา รัตนสุภา)
ผู้อำนวยการ
ศูนย์วิจัยและเตือนภัยมลพิษโรงงานภาคใต้