

ภาคผนวก ง
หนังสือรับรองการขึ้นทะเบียน
ห้องปฏิบัติการวิเคราะห์เอกชน

ที่ อก ๐๓๑๐(๑)/ ๕๖ ๙ ๑



กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๕๐๐

๐๗ กรกฎาคม ๒๕๖๕

เรื่อง เปลี่ยนแปลงบุคลากรและสารมลพิษที่วิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท ยูนิเทค แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๙ พฤษภาคม ๒๕๖๔

สิ่งที่ส่งมาด้วย เอกสารแนบท้ายหนังสือเปลี่ยนแปลงบุคลากรและสารมลพิษที่วิเคราะห์
บริษัท ยูนิเทค แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด จำนวน ๒ แผ่น

ตามหนังสือที่อ้างถึง บริษัท ยูนิเทค แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด
ห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ๖-๑๔๕ สถานที่ตั้งเลขที่ ๓ ซอยอุดมสุข ๔๑ ถนนสุขุมวิท
แขวงบางจาก เขตพระโขนง กรุงเทพมหานคร ขอเปลี่ยนแปลงบุคลากรและสารมลพิษที่วิเคราะห์ที่
ความละเอียดแจ้งแล้ว นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว มีความเห็นดังนี้

๑. ให้อยกเลิกเจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกชน จำนวน ๓ ราย

๒. ให้เพิ่มผู้ควบคุมห้องปฏิบัติการวิเคราะห์เอกชน จำนวน ๑ ราย

๓. ให้เพิ่มขอบข่ายสารมลพิษที่วิเคราะห์ในภาคดิน อากาศเสีย และดิน ตามสิ่งที่ส่งมาด้วย

อนึ่ง หนังสือฉบับนี้จะหมดอายุพร้อมหนังสือขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

ในวันที่ ๒ กุมภาพันธ์ ๒๕๗๒

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

๖๐๘ ๙๙

(นายประสม คำวังพวง)

ผู้อำนวยการกองวิจัยและเฝ้าระวังมลพิษโรงงาน
ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม

กองวิจัยและเฝ้าระวังมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕

โทรสาร ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๔๙

ไปรษณีย์อิเล็กทรอนิกส์ saraban@diw.mail.go.th



อำนาจถูกต้อง



"อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว"



- ๒ -

3. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Acid Digestion of Sediments, Sludges, and Soils. SW-846 Method 3050B, 1996.

4. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Volatile Organic Compounds in Various Sample Matrices Using Equilibrium Headspace Analysis. SW-846 Method 5021A, 2014.

5. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Inductively Coupled Plasma-Optical Emission Spectrometry. SW-846 Method 6010D, 2018.

6. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Flame Atomic Absorption Spectrophotometry. SW-846 Method 7000B, 2007.

7. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Volatile Organic Compounds by Gas Chromatography/Mass Spectrometry. SW-846 Method 8260D, 2018.

8. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Soil and Waste pH. SW-846 Method 9045D, 2004.

อำนาจ

เอกสารแนบท้ายหนังสือเปลี่ยนแปลงบุคลากรและสารมลพิษที่วิเคราะห์

บริษัท ยูนิเทค แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด เลขทะเบียน ๖-๑๔๕

ที่ อก ๐๓๑๐(๑)/ ๕๖ ๙ ๑

ลงวันที่ ๐๗ กรกฎาคม ๒๕๖๕

ขอขยายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๑๑ รายการ

น้ำได้ดิน จำนวน 4 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Aluminum	Digestion, Inductively Coupled Plasma Method ^[1]
2	Copper	1) Digestion, Flame Atomic Absorption Spectrometric Method ^[1] 2) Digestion, Inductively Coupled Plasma Method ^[1]
3	Iron	1) Digestion, Flame Atomic Absorption Spectrometric Method ^[1] 2) Digestion, Inductively Coupled Plasma Method ^[1]
4	Molybdenum	Digestion, Inductively Coupled Plasma Method ^[1]

อากาศเสีย (ปล่องระบาย) จำนวน 1 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Oxides of Nitrogen	Absorption Sampling, Ion Chromatographic Method ^[2]

ดิน จำนวน 6 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Aluminum	Digestion, Inductively Coupled Plasma Method ^[3,5]
2	Copper	1) Digestion, Flame Atomic Absorption Spectrometric Method ^[3,6] 2) Digestion, Inductively Coupled Plasma Method ^[3,5]
3	Iron	1) Digestion, Flame Atomic Absorption Spectrometric Method ^[3,6] 2) Digestion, Inductively Coupled Plasma Method ^[3,5]
4	Molybdenum	Digestion, Inductively Coupled Plasma Method ^[3,6]
5	pH	Electrometric Method ^[8]
6	TPH (C ₅ -C ₆)	Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^[4,7]

เอกสารอ้างอิง

1. APHA, AWWA, WEF. Standard Methods for the Examination of Water and Wastewater. 24th ed. Washington, DC: APHA, 2023.

2. United States Environmental Protection Agency. Standards of Performance for New Stationary Sources. 40 CFR 60. Appendix A, 2023.

3. United States...



ที่ อก ๐๓๑๐(๑)/ ๑๔ ๙ ๑

กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๕๐๐

๒๔ กุมภาพันธ์ ๒๕๖๕

เรื่อง ยกเลิกบุคลากรของห้องปฏิบัติการวิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท ยูนิเทค แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๑๘ กุมภาพันธ์ ๒๕๖๔

ตามคำขอที่อ้างถึง บริษัท ยูนิเทค แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด
ห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ๖-๑๔๕ สถานที่ตั้งเลขที่ ๓ ซอยอุดมสุข ๔๑ ถนนสุขุมวิท
แขวงบางจาก เขตพระโขนง กรุงเทพมหานคร ขอยกเลิกบุคลากร ความละเอียดแจ้งแล้ว นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้อยกเลิกเจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกชน
จำนวน ๓ ราย ได้แก่



จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

รังสรรค์

(นายธีรศักดิ์ อิศรางกูร ณ อยุธยา)

รองอธิบดี ปฏิบัติราชการแทน

อธิบดีกรมโรงงานอุตสาหกรรม

กองวิจัยและเฝ้าระวังมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕

โทรสาร ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๔๙

ไปรษณีย์อิเล็กทรอนิกส์ saraban@diw.mail.go.th



อำนาจถูกต้อง



ที่ อก ๐๓๑๐(๑)/ ๑ ๐ ๘ ๘



กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๕๐๐

๐๗ กุมภาพันธ์ ๒๕๖๘

เรื่อง ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

เรียน กรรมการผู้จัดการ บริษัท ยูโนเด็ค แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๓ ธันวาคม ๒๕๖๗

- สิ่งที่ส่งมาด้วย ๑. รายชื่อผู้ควบคุมห้องปฏิบัติการวิเคราะห์เอกชน จำนวน ๔๐ ราย
๒. รายชื่อเจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกชน จำนวน ๑๔๑ ราย
๓. ขอบข่ายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม

ตามคำขอที่อ้างถึง บริษัท ยูโนเด็ค แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด ขอต่ออายุ
หนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ๖-๑๔๕ สถานที่ตั้งเลขที่ ๓ ซอยอุดมสุข ๔๑
ถนนสุขุมวิท แขวงบางจาก เขตพระโขนง กรุงเทพมหานคร ต่อกรมโรงงานอุตสาหกรรม นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้บริษัท ยูโนเด็ค แอนนาลิสต์ แอนด์ เอ็นจิเนียริง
คอนซัลแตนท์ จำกัด ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน โดยมีองค์ประกอบดังนี้

- ก. ผู้ควบคุมห้องปฏิบัติการวิเคราะห์เอกชน จำนวน ๔๐ ราย ตามสิ่งที่ส่งมาด้วย ๑
ข. เจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกชน จำนวน ๑๔๑ ราย ตามสิ่งที่ส่งมาด้วย ๒
ค. ขอบข่ายสารมลพิษที่ได้รับขึ้นทะเบียนให้วิเคราะห์ในน้ำ/น้ำเสีย น้ำใต้ดิน อากาศเสีย
สิ่งปฏิกูลหรือวัสดุที่ไม่ใช้แล้ว และดิน ตามสิ่งที่ส่งมาด้วย ๓

หนังสือฉบับนี้จะหมดอายุในวันที่ ๒ กุมภาพันธ์ ๒๕๗๒ หากประสงค์จะต่ออายุหนังสือ
รับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ให้ยื่นคำขอต่ออายุพร้อมเอกสารประกอบคำขอต่อกรมโรงงาน
อุตสาหกรรมภายใน ๖๐ วัน ก่อนวันสิ้นสุดอายุของหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

อธิบดีกรม
(นายธีรศักดิ์ อิศรางกูร ณ อยุธยา)

รองอธิบดี ปฏิบัติราชการแทน

อธิบดีกรมโรงงานอุตสาหกรรม

UNITED ANALYST AND ENGINEERING
CONSULTANT COMPANY LIMITED

ดำเนินการถูกต้อง

กองวิจัยและเตือนภัยมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๔๓๐ ๖๓๓๒ ต่อ ๒๑๐๓-๕

โทรสาร ๐ ๒๔๓๐ ๖๓๓๒ ต่อ ๒๑๕๔

ไปรษณีย์อิเล็กทรอนิกส์ sarabandiw@mail.go.th



"อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว"



เอกสารแนบท้ายหนังสือต่ออายุรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท ยูโนเด็ค แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด เลขทะเบียน ๖-๑๔๕
ที่ อก ๐๓๑๐(๑)/ ๑ ๐ ๘ ๘ ลงวันที่ ๐๗ กุมภาพันธ์ ๒๕๖๘

ก. ผู้ควบคุมห้องปฏิบัติการวิเคราะห์เอกชน จำนวน ๔๐ ราย



ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๐๑
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๐๒
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๐๓
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๐๔
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๐๕
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๐๖
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๐๗
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๐๘
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๐๙
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๑๐
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๑๑
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๑๒
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๑๓
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๑๔
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๑๕
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๑๖
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๑๗
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๑๘
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๑๙
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๒๐
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๒๑
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๒๒
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๒๓
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๒๔
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๒๕
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๒๖
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๒๗
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๒๘
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๒๙
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๓๐
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๓๑
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๓๒
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๓๓
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๓๔
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๓๕
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๓๖
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๓๗
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๓๘
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๓๙
ทะเบียนเลขที่ ๖-๑๔๕-ก-๐๐๔๐


ตามนี้

๓๖) นายนาเคนทร์...

เอกสารแนบท้ายหนังสือต่ออายุรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท ยูโนเด็ค แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด เลขทะเบียน ๖-๑๔๕
ที่ อก ๐๓๑๐(๑)/ ๑ ๐ ๘ ๘ ลงวันที่ ๐๗ กุมภาพันธ์ ๒๕๖๘

ข. เจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกชน จำนวน ๑๔๑ ราย

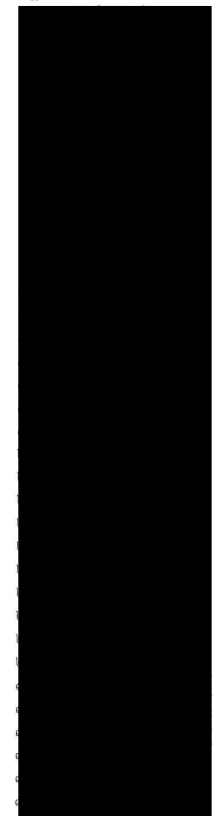

ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๐๑
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๐๒
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๐๓
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๐๔
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๐๕
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๐๖
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๐๗
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๐๘
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๐๙
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๑๐
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๑๑
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๑๒
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๑๓
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๑๔
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๑๕
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๑๖
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๑๗
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๑๘
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๑๙
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๒๐
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๒๑
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๒๒
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๒๓
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๒๔
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๒๕
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๒๖
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๒๗
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๒๘
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๒๙
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๓๐
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๓๑
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๓๒
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๓๓
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๓๔
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๓๕
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๓๖
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๓๗
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๓๘
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๓๙
ทะเบียนเลขที่ ๖-๑๔๕-ข-๐๐๔๐

ตามนี้

๓๖) นางสาวนิภาพร...

UAE
UNITED ANALYST AND ENGINEERING
CONSULTANT COMPANY LIMITED

ดำเนินการถูกต้อง

UAE
UNITED ANALYST AND ENGINEERING
CONSULTANT COMPANY LIMITED

9252

ทะเบียนเลขที่ 7-๑๔๕-จ-๐๑๔๐ มงชภัท
ทะเบียนเลขที่ 7-๑๔๕-จ-๐๑๔๑
ทะเบียนเลขที่ 7-๑๔๕-จ-๐๑๔๒

၇၁၂

ลำดับ	สารเคมี	วิธีวิเคราะห์
1	Aldrin	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
2	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
3	Barium	Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
4	α -BHC	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
5	β -BHC	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
6	δ -BHC	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
7	γ -BHC	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
8	Biochemical Oxygen Demand	1) 5-Day BOD Test, Azide Modification Method ⁽⁴⁾ 2) 5-Day BOD Test, Membrane Electrode Method ⁽⁴⁾
9	Cadmium	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
10	Chemical Oxygen Demand	1) Closed Reflux, Titrimetric Method ⁽⁴⁾ 2) Closed Reflux, Colorimetric Method ⁽⁴⁾ 3) Open Reflux, Titrimetric Method ⁽⁴⁾
11	Chlordane	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
12	Chromium	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
13	Color	ADMI Weighted-Ordinate Spectrophotometric Method ⁽⁴⁾
14	Copper	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
15	Cyanide	1) Distillation, Colorimetric Method ⁽⁴⁾ 2) Total Cyanide after Distillation, by Flow Injection Analysis Method ⁽⁴⁾
16	<i>o,p'</i> -DDT	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
17	4,4'-DDD	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
18	4,4'-DDE	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
19	4,4'-DDT	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
20	Dieldrin	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
21	Endosulfan I	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
22	Endosulfan II	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
23	Endosulfan sulfate	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
24	Endrin	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾

25 Endrin aldehyde...

UAE
UNITED ANALYST AND ENGINEERING
CONSULTANT COMPANY LIMITED

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
25	Endrin aldehyde	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
26	Formaldehyde	Distillation, Colorimetric Method ⁽²⁾
27	Free Chlorine	1) Iodometric Method ⁽⁴⁾ 2) DPD Ferrous Titrimetric Method ⁽⁴⁾
28	Heptachlor	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
29	Heptachlor Epoxide	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
30	Hexavalent Chromium	Colorimetric Method ⁽⁴⁾
31	Lead	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
32	Manganese	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
33	Mercury	Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ⁽⁴⁾
34	Methoxychlor	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
35	Nickel	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
36	Oil & Grease	1) Liquid-Liquid, Partition-Gravimetric Method ⁽⁴⁾ 2) Soxhlet Extraction Method ⁽⁴⁾
37	pH	Electrometric Method ⁽⁴⁾
38	Phenols	1) Distillation, Chloroform Extraction Method ⁽⁴⁾ 2) Distillation, Direct Photometric Method ⁽⁴⁾
39	Selenium	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
40	Sulfide	1) Iodometric Method ⁽⁴⁾ 2) Methylene Blue Method ⁽⁴⁾
41	Temperature	Laboratory and Field Methods ⁽⁴⁾
42	Total Dissolved Solids	Dried at 180 °C ⁽⁴⁾
43	Total Kjeldahl Nitrogen	Semi-Micro-Kjeldahl Method ⁽⁴⁾
44	Total Suspended Solids	Dried from 103 to 105 °C ⁽⁴⁾
45	Trivalent Chromium	1) Digestion, Direct Air-Acetylene Flame Method; ^{วงษ์ภัทรา} Colorimetric Method; Calculation ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method; ^{วงษ์ภัทรา} Colorimetric Method; Calculation ⁽⁴⁾
46	Zinc	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾

น้ำได้คืน...

น้ำได้คืน จำนวน 126 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Acenaphthene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
2	Acetone	2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾ Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
3	Aldrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
4	Anthracene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
5	Antimony	Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
6	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
7	Atrazine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
8	Barium	Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
9	Benz(a)anthracene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
10	Benzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
11	Benzo(b)fluoranthene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
12	Benzo(k)fluoranthene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
13	Benzoic acid	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

14 Benzo(a)pyrene...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
14	Benzo(a)pyrene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
15	Benzo(g,h,i)perylene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
16	Beryllium	Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
17	Bis(2-chloroethyl)ether	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
18	Bis(2-ethylhexyl)phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
19	Bromodichloromethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
20	Bromoform	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
21	Butanol	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
22	Butyl benzyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
23	Cadmium	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ⁽⁴⁾ 3) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
24	Carbazole	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
25	Carbon disulfide	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
26	Carbon tetrachloride	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
27	Chlordane	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
28	p-Chloroaniline	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

29 Chlorobenzene...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
29	Chlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
30	Chlorodibromomethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
31	Chloroform	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
32	2-Chlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
33	Chromium	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
34	Chromium (III)	1) Digestion, Direct Air-Acetylene Flame Method; Colorimetric Method; Calculation ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method; Colorimetric Method; Calculation ⁽⁴⁾
35	Chromium (VI)	Colorimetric Method ⁽⁴⁾
36	Chrysene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
37	Cyanide	Distillation, Colorimetric Method ⁽⁴⁾
38	2,4-D	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
39	DDD	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
40	DDE	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
41	DDT	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
42	Dibenz(a,h)anthracene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

43 Di-n-butyl phthalate...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
43	Di-n-butyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
44	1,2-Dichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
45	1,3-Dichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
46	1,4-Dichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
47	3,3'-Dichlorobenzidine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
48	1,1-Dichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
49	1,2-Dichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
50	1,1-Dichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
51	cis-1,2-Dichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
52	trans-1,2-Dichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
53	2,4-Dichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
54	1,2-Dichloropropane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
55	1,3-Dichloropropane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
56	1,3-Dichloropropene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
57	Dieldrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
58	Diethyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
59	2,4-Dimethylphenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
60	2,4-Dinitrophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

61 2,4-Dinitrotoluene...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
61	2,4-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
62	2,6-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
63	Di-n-Octyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
64	Endosulfan	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
65	Endrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
66	Ethylbenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
67	Fluoranthene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
68	Fluorene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
69	Heptachlor	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
70	Heptachlor epoxide	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
71	Hexachlorobenzene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
72	Hexachloro-1,3-butadiene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
73	n-Hexane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

74 α -HCH...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
74	α -HCH	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
75	β -HCH	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
76	γ -HCH	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
77	Hexachlorocyclopentadiene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
78	Hexachloroethane	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
79	Indeno(1,2,3-cd)pyrene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
80	Isophorone	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
81	Lead	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ⁽⁴⁾ 3) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
82	Manganese	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
83	Mercury	Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ⁽⁴⁾
84	Methanol	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
85	Methoxychlor	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
86	Methyl bromide	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

87 Methylene chloride...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
87	Methylene chloride	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
88	2-Methylphenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
89	2-Methylnaphthalene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
90	Methyl tert-butyl ether	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
91	Naphthalene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
92	Nickel	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
93	Nitrobenzene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
94	N-Nitrosodiphenylamine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
95	N-Nitrosodi-n-propylamine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
96	Polychlorinated Biphenyls - PCB 1016 - PCB 1221 - PCB 1232 - PCB-1242 - PCB-1248 - PCB-1254 - PCB-1260	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
97	Pentachlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
98	pH	Electrometric Method ⁽⁴⁾
99	Phenanthrene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

100 Phenol...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
100	Phenol	1) Distillation, Chloroform Extraction Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
101	Pyrene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
102	Selenium	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
103	Silver	Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
104	Styrene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
105	1,1,2,2-Tetrachloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
106	Tetrachloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
107	Toluene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
108	Toxaphene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
109	TPH (C ₅ - C ₆)	1) Purge and Trap, Gas Chromatographic Method ^(12,22) 2) Purge and Trap, Gas Chromatographic/Mass spectrometric Method ^(12,27)
110	TPH (C ₈ - C ₁₆)	Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(9,22)
111	TPH (C ₁₆ - C ₃₅)	Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(9,22)
112	1,2,4-Trichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
113	1,1,1-Trichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
114	1,1,2-Trichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
115	Trichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

116 2,4,5-Trichlorophenol...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
116	2,4,5-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
117	2,4,6-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
118	1,3,5-Trimethylbenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
119	Vanadium	Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
120	Vinyl acetate	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
121	Vinyl chloride	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
122	m-Xylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
123	o-Xylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
124	p-Xylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
125	Xylene (Total)	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
126	Zinc	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾

อากาศเสีย (ปล่องระบาย) จำนวน 25 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Antimony	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽⁵⁾
2	Arsenic	1) Isokinetic Sampling, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ⁽⁵⁾ 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽⁵⁾
3	Cadmium	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ⁽⁵⁾ 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽⁵⁾
4	Carbon Monoxide	Instrumental Analyzer Method ⁽⁵⁾
5	Chlorine	Isokinetic Sampling, Ion Chromatographic Method ⁽⁵⁾
6	Chromium	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ⁽⁵⁾

Chromium (ต่อ)...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
6	Chromium (ต่อ)	2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽⁵⁾
7	Cobalt	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽⁵⁾
8	Copper	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ⁽⁵⁾ 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽⁵⁾
9	Cresol	Absorption Sampling, Gas Chromatographic Method ⁽⁵⁾
10	Dioxins/Furans	Isokinetic Sampling ⁽⁵⁾
11	Hydrogen Chloride	Isokinetic Sampling, Ion Chromatographic Method ⁽⁵⁾
12	Hydrogen Fluoride	Isokinetic Sampling, Ion Chromatographic Method ⁽⁵⁾
13	Hydrogen Sulfide	Absorption Sampling, Iodometric Method ⁽⁵⁾
14	Lead	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ⁽⁵⁾ 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽⁵⁾
15	Manganese	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ⁽⁵⁾ 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽⁵⁾
16	Mercury	Isokinetic Sampling, Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ⁽⁵⁾
17	Nickel	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ⁽⁵⁾ 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽⁵⁾
18	Opacity	Ringelmann's Method ⁽¹⁾
19	Oxides of Nitrogen	1) Absorption Sampling, Phenoldisulfonic acid Method ⁽⁵⁾ 2) Instrumental Analyzer Method ⁽⁵⁾
20	Selenium	1) Isokinetic Sampling, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ⁽⁵⁾ 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽⁵⁾
21	Sulfur Dioxide	1) Absorption Sampling, Barium-Thorin Titrimetric Method ⁽⁵⁾ 2) Instrumental Analyzer Method ⁽⁵⁾
22	Sulfuric Acid	Isokinetic Sampling, Barium-Thorin Titrimetric Method ⁽⁵⁾

23 Total Suspended Particulate...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
23	Total Suspended Particulate	Isokinetic Sampling, Gravimetric Method ⁽⁵⁾
24	Vanadium	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽⁵⁾
25	Xylene	1) Bag Sampling, Gas Chromatographic Method ⁽⁵⁾ 2) Adsorption Sampling, Gas Chromatographic Method ⁽⁵⁾

สิ่งปฏิกูลหรือวัสดุที่ไม่ใช้แล้ว จำนวน 35 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Aldrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
2	Antimony	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
3	Arsenic	1) Waste Extraction, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^(3,6,16) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 3) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^(7,14) 4) Digestion, Inductively Coupled Plasma Method ^(7,14)
4	Barium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
5	Beryllium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
6	Cadmium	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^(3,6,15) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 3) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,14) 4) Digestion, Inductively Coupled Plasma Method ^(7,14)
7	Chlordane	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)

8 Chromium...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
8	Chromium	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^(3,6,15) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 3) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 4) Digestion, Inductively Coupled Plasma Method ^(7,14)
9	Chromium (III)	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method; Waste Extraction, Colorimetric Method; Calculation ^(3,6,15,17) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method; Waste Extraction, Colorimetric Method; Calculation ^(3,6,14,17) 3) Digestion, Flame Atomic Absorption Spectrometric Method; Alkaline Digestion, Colorimetric Method; Calculation ^(7,8,15,17) 4) Digestion, Inductively Coupled Plasma Method; Alkaline Digestion, Colorimetric Method; Calculation ^(7,8,14,17)
10	Chromium (VI)	1) Waste Extraction, Colorimetric Method ^(3,17) 2) Alkaline Digestion, Colorimetric Method ^(8,17)
11	Cobalt	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
12	Copper	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^(3,6,15) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 3) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 4) Digestion, Inductively Coupled Plasma Method ^(7,14)
13	2,4-D	1) Waste Extraction, Gas Chromatographic Method ^(3,26) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
14	DDD	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)

15 DDE...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
15	DDE	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
16	DDT	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
17	Dieldrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
18	Endrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
19	Heptachlor	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
20	Lead	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^(3,6,15) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 3) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 4) Digestion, Inductively Coupled Plasma Method ^(7,14)
21	Lindane	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
22	Mercury	1) Waste Extraction, Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^(3,19) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 3) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ⁽¹⁹⁾ 4) Digestion, Inductively Coupled Plasma Method ^(7,14)

Mercury (ต่อไป)...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
22	Mercury (ต่อไป)	5) Thermal Decomposition Amalgamation and Atomic Absorption Spectrometric Method ⁽²⁰⁾
23	Methoxychlor	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
24	Molybdenum	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
25	Nickel	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^(3,6,15) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 3) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 4) Digestion, Inductively Coupled Plasma Method ^(7,14)
26	Polychlorinated Biphenyls - Aroclor 1016 - Aroclor 1221 - Aroclor 1232 - Aroclor 1242 - Aroclor 1248 - Aroclor 1254 - Aroclor 1260 - 2-Chlorobiphenyl - 2,3-Dichlorobiphenyl - 2,2',5'-Trichlorobiphenyl - 2,4',5'-Trichlorobiphenyl - 2,2',3,5'-Tetrachlorobiphenyl - 2,2',5,5'-Tetrachlorobiphenyl - 2,3',4,4'-Tetrachlorobiphenyl - 2,2',3,4,5'-Pentachlorobiphenyl - 2,2',4,5,5'-Pentachlorobiphenyl - 2,3,3',4',6'-Pentachlorobiphenyl	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,24) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,24)

Polychlorinated Biphenyls(ต่อไป)...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
	Polychlorinated Biphenyls(ต่อไป) - 2,2',3,4,4',5'-Hexachlorobiphenyl - 2,2',3,4,5,5'-Hexachlorobiphenyl - 2,2',3,5,5',6'-Hexachlorobiphenyl - 2,2',4,4',5,5'-Hexachlorobiphenyl - 2,2',3,3',4,4',5'-Heptachlorobiphenyl - 2,2',3,4,4',5,5'-Heptachlorobiphenyl - 2,2',3,4,4',5',6'-Heptachlorobiphenyl - 2,2',3,4',5,5',6'-Heptachlorobiphenyl - 2,2',3,3',4,4',5,5',6'-Nonachlorobiphenyl	
27	Pentachlorophenol	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^(3,9,28) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
28	pH	Electrometric Method ^(31,32)
29	Selenium	1) Waste Extraction, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^(3,6,21) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 3) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^(7,21) 4) Digestion, Inductively Coupled Plasma Method ^(7,14)
30	Silver	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
31	Thallium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)

32 Toxaphene...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
32	Toxaphene	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
33	Trichloroethylene	1) Waste Extraction, Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(3,12,27) 2) Waste Extraction, Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(3,11,27) 3) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 4) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
34	Vanadium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
35	Zinc	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^(3,6,15) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 3) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 4) Digestion, Inductively Coupled Plasma Method ^(7,14)

ดิน จำนวน 125 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Acenaphthene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
2	Acetone	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
3	Aldrin	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
4	Anthracene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25)

Anthracene (ต่อ)...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
4	Anthracene (ต่อ)	2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
5	Antimony	Digestion, Inductively Coupled Plasma Method ^(7,14)
6	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^(7,16) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
7	Atrazine	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
8	Barium	Digestion, Inductively Coupled Plasma Method ^(7,14)
9	Benz(a)anthracene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
10	Benzene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
11	Benzo(b)fluoranthene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
12	Benzo(k)fluoranthene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
13	Benzoic acid	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
14	Benzo(a)pyrene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
15	Benzo(g,h,i)perylene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
16	Beryllium	Digestion, Inductively Coupled Plasma Method ^(7,14)

17 Bis(2-chloroethyl)ether...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
17	Bis(2-chloroethyl)ether	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
18	Bis(2-ethylhexyl)phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
19	Bromodichloromethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
20	Bromoform	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
21	Butanol	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
22	Butyl benzyl phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
23	Cadmium	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
24	Carbazole	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
25	Carbon disulfide	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
26	Carbon tetrachloride	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
27	Chlordane	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
28	p-Chloroaniline	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
29	Chlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
30	Chlorodibromomethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
31	Chloroform	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
32	2-Chlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)

33 Chromium...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
33	Chromium	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
34	Chromium (III)	1) Digestion, Flame Atomic Absorption Spectrometric Method; Alkaline Digestion, Colorimetric Method; Calculation ^(7,8,15,17) 2) Digestion, Inductively Coupled Plasma Method; Alkaline Digestion, Colorimetric Method; Calculation ^(7,8,14,17)
35	Chromium (VI)	Alkaline Digestion, Colorimetric Method ^(8,17)
36	Chrysene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
37	Cyanide	Extraction, Distillation, Colorimetric Method ^(29,30)
38	2,4-D	Ultrasonic Extraction, Gas Chromatographic Method ⁽²⁶⁾
39	DDD	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
40	DDE	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
41	DDT	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
42	Dibenz(a,h)anthracene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
43	Di-n-butyl phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
44	1,2-Dichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)

45 1,3-Dichlorobenzene...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
45	1,3-Dichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
46	1,4-Dichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
47	3,3'-Dichlorobenzidine	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
48	1,1-Dichloroethane	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
49	1,2-Dichloroethane	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
50	1,1-Dichloroethylene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
51	cis-1,2-Dichloroethylene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
52	trans-1,2-Dichloroethylene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
53	2,4-Dichlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
54	1,2-Dichloropropane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
55	1,3-Dichloropropane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
56	1,3-Dichloropropene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
57	Dieldrin	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)

58 Diethyl phthalate...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
58	Diethyl phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
59	2,4-Dimethylphenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
60	2,4-Dinitrophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
61	2,4-Dinitrotoluene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
62	2,6-Dinitrotoluene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
63	Di-n-Octyl phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
64	Endosulfan	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
65	Endrin	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
66	Ethylbenzene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
67	Fluoranthene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
68	Fluorene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
69	Heptachlor	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
70	Heptachlor epoxide	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)


Heptachlor epoxide (ต่อ)...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
70	Heptachlor epoxide (ต่อ)	2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
71	Hexachlorobenzene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
72	Hexachloro-1,3-butadiene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
73	n-Hexane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
74	α-HCH	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
75	β-HCH	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
76	γ-HCH	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
77	Hexachlorocyclopentadiene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
78	Hexachloroethane	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
79	Indeno(1,2,3-cd)pyrene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
80	Isophorone	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
81	Lead	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
82	Manganese	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)

83 Mercury...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
83	Mercury	1) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ⁽¹⁹⁾ 2) Thermal Decomposition Amalgamation and Atomic Absorption Spectrometric Method ⁽²⁰⁾
84	Methanol	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
85	Methoxychlor	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
86	Methyl bromide	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
87	Methylene chloride	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
88	2-Methylphenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
89	2-Methylnaphthalene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
90	Methyl tert-butyl ether	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
91	Naphthalene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
92	Nickel	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
93	Nitrobenzene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
94	N-Nitrosodiphenylamine	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
95	N-Nitrosodi-n-propylamine	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
96	Polychlorinated Biphenyls - Aroclor 1016	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,24)

Polychlorinated Biphenyls(ต่อ)...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
96	Polychlorinated Biphenyls(ต่อ) - Aroclor 1221 - Aroclor 1232 - Aroclor 1242 - Aroclor 1248 - Aroclor 1254 - Aroclor 1260 Polychlorinated Biphenyls - 2-Chlorobiphenyl - 2,3-Dichlorobiphenyl - 2,2',5-Trichlorobiphenyl - 2,4',5-Trichlorobiphenyl - 2,2',3,5'-Tetrachlorobiphenyl - 2,2',5,5'-Tetrachlorobiphenyl - 2,3',4,4'-Tetrachlorobiphenyl - 2,2',3,4,5'- Pentachlorobiphenyl - 2,2',4,5,5'- Pentachlorobiphenyl - 2,3,3',4',6- Pentachlorobiphenyl - 2,2',3,4,4',5'- Hexachlorobiphenyl - 2,2',3,4,5,5'- Hexachlorobiphenyl - 2,2',3,5,5',6- Hexachlorobiphenyl - 2,2',4,4',5,5'- Hexachlorobiphenyl - 2,2',3,3',4,4',5- Heptachlorobiphenyl - 2,2',3,4,4',5,5'- Heptachlorobiphenyl - 2,2',3,4,4',5,6- Heptachlorobiphenyl - 2,2',3,4',5,5',6- Heptachlorobiphenyl - 2,2',3,3',4,4',5,5',6- Nonachlorobiphenyl	2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28) Ultrasonic Extraction, Gas Chromatographic Method ^(10,24)  ดำเนินการถูกต้อง อนุมัติ

97 Pentachlorophenol...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
97	Pentachlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
98	Phenanthrene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
99	Phenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
100	Pyrene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
101	Selenium	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^(7,21) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
102	Silver	Digestion, Inductively Coupled Plasma Method ^(7,14)
103	Styrene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
104	1,1,2,2-Tetrachloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
105	Tetrachloroethylene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
106	Toluene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
107	Toxaphene	Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
108	TPH (C ₅ -C ₈)	1) Purge and Trap, Gas Chromatographic Method ^(13,27) 2) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
109	TPH (C ₈ -C ₁₆)	Ultrasonic Extraction, Gas Chromatographic Method ^(10,22)
110	TPH (C ₁₆ -C ₃₅)	Ultrasonic Extraction, Gas Chromatographic Method ^(10,22)

111 1,2,4-Trichlorobenzene...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
111	1,2,4-Trichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
112	1,1,1-Trichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
113	1,1,2-Trichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
114	Trichloroethylene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
115	2,4,5-Trichlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
116	2,4,6-Trichlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
117	1,3,5-Trimethylbenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
118	Vanadium	Digestion, Inductively Coupled Plasma Method ^(7,14)
119	Vinyl acetate	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
120	Vinyl chloride	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
121	m-Xylene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
122	o-Xylene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
123	p-Xylene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
124	Xylene (Total)	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)

125 Zinc...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
125	Zinc	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,13) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)

เอกสารอ้างอิง

- กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2549. เรื่อง กำหนดค่าปริมาณเขม่าควันที่เจือปนในอากาศที่ระบายออกจากปล่องของหม้อน้ำโรงสีข้าวที่ใช้แก๊สเป็นเชื้อเพลิง. ราชกิจจานุเบกษา. 4 ธันวาคม 2549. เล่มที่ 123 ตอนพิเศษ 125 ง.
- สมาคมวิศวกรรมสิ่งแวดล้อมแห่งประเทศไทย. คู่มือวิเคราะห์น้ำเสีย. พิมพ์ครั้งที่ 4. กรุงเทพฯ: เรือนแก้วกรพิมพ์. 2547.
- กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2566. เรื่อง การจัดการสิ่งปฏิกูลหรือวัสดุที่ไม่ใช้แล้ว. ราชกิจจานุเบกษา. 31 พฤษภาคม 2566. เล่มที่ 140 ตอนพิเศษ 126 ง.
- APHA, AWWA, WEF. Standard Methods for the Examination of Water and Wastewater. 24th ed. Washington, DC: APHA, 2023.
- United States Environmental Protection Agency. Standards of Performance for New Stationary Sources. 40 CFR 60. Appendix A, 2020.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. SW-846, 2014.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Acid Digestion of Sediments, Sludges, and Soils. SW-846 Method 3050B, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Alkaline Digestion for Hexavalent Chromium. SW-846 Method 3060A, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste 3. Physical/Chemical Methods. Separatory Funnel Liquid-Liquid Extraction. SW-846 Method 3510C, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Ultrasonic Extraction. SW-846 Method 3550C, 2007.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Volatile Organic Compounds in Various Sample Matrices Using Equilibrium Headspace Analysis. SW-846 Method 5021A, 2014.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Purge and Trap for Aqueous Samples. SW-846 Method 5030C, 2003.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Closed System Purge and Trap and Extraction for Volatile Organics in Soil and Waste Sample. SW-846 Method 5035A, 2000.
- United States...

14. United States...

14. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Inductively Coupled Plasma-Optical Emission Spectrometry. SW-846 Method 6010D, 2014.

15. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Flame Atomic Absorption Spectrophotometry. SW-846 Method 7000B, 2007.

16. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Arsenic (Atomic Absorption, Gaseous Hydride). SW-846 Method 7061A, 1992.

17. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Chromium, Hexavalent (Colorimetric). SW-846 Method 7196A, 1992.

18. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Mercury in Liquid Waste (Manual Cold Vapor Technique). SW-846 Method 7470A, 1994.

19. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Mercury in Solid or Semisolid Waste (Manual Cold-Vapor Technique). SW-846 Method 7471B, 1998.

20. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Mercury in Solids and Solutions by Thermal Decomposition, Amalgamation, and Atomic Absorption Spectrophotometry. SW-846 Method 7473, 2007.

21. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Selenium (Atomic Absorption, Borohydride Reduction). SW-846 Method 7742, 1994.

22. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Nonhalogenated Organics Using GC/FID. SW-846 Method 8015D, 2003.

23. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Organochlorine Pesticides by Gas Chromatography. SW-846 Method 8081B, 2007.

24. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Polychlorinated Biphenyls (PCBs) by Gas Chromatography. SW-846 Method 8082A, 2007.

25. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Organochlorine Pesticides by Gas Chromatography. SW-846 Method 8081B, 2007.

26. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Chlorinated Herbicides by GC Using Methylation or Pentafluorobenzoylation Derivatization. SW-846 Method 8151A, 1996.

27. United States...

27. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Volatile Organic Compounds by Gas Chromatography/ Mass Spectrometry. SW-846 Method 8260D, 2018.

28. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Semivolatile Organic Compounds by Gas Chromatography/Mass Spectrometry. SW-846 Method 8270E, 2018.

29. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Cyanide Extraction Procedure for Solids and Oils. SW-846 Method 9013A, 2014.

30. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Cyanide in Waters and Extracts using Titrimetric and Manual Spectrophotometric Procedures. SW-846 Method 9014, 2014.

31. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. pH Electrometric Measurement. SW-846 Method 9040C, 2004.

32. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Soil and Waste pH. SW-846 Method 9045D, 2004.

๓๑๑



ภาคผนวก จ
เอกสารสอบเทียบเครื่องมือ

List of Instruments Certification for Air Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration	Remark
Ambient									
1	Standard Gases (Mixture)	Nitrogen Dioxide Sulphur Dioxide Carbon Monoxide	Airgas	EB0159156 2015PSIG	Airgas an Air Liquide company	E04NJ99E15A01D3	6 Nov 23	6 Nov 26	-
2	Standard Gas	Total Hydrocarbons	Air Liquide	CC143232	Airgas an Air Liquide company	E03A99E15A006C	16 Oct 20	16 Oct 28	-
3	Nitrogen Dioxide Analyzer	Nitrogen Dioxide	Thermo Scientific	42i CM22387040	UAE Consultant Co.,Ltd.	20092024	20 Sep 24	19 Sep 25	-
5	Sulphur Dioxide Analyzer	Sulphur Dioxide	Thermo Scientific	43i 1180540065	UAE Consultant Co.,Ltd.	4092024	4 Sep 24	4 Sep 25	-
6	Carbon Monoxide Analyzer	Carbon Monoxide	Thermo Scientific	48i CM08140004	UAE Consultant Co.,Ltd.	03092024	3 Sep 24	3 Sep 25	-
7	Total Hydrocarbons Analyzer	Total Hydrocarbons	Thermo Scientific	APHA-370 93JN1MN9	UAE Consultant Co.,Ltd.	2092024	2 Sep 24	2 Sep 25	-
Water									
8	pH Merer	pH Meter	HANNA Instruments	HI2020-02 C0051107	National Food Institute, Ministry of Industry, Thailand	2404042-001-01	14 Aug 24	13 Aug 25	-
9	Balance	ของแข็งแขวนลอยทั้งหมด (Total Suspended Solids) ของแข็งละลายน้ำทั้งหมด (Total Dissolved Solids) น้ำมันและไขมัน (Oil & Grease) ตะกอนหนัก	Mettler-Toledo	XPE205/ B748058497	Mettler-Toledo (Thailand) Ltd.	TH2084-013-092024-ACC-TH	21 Sep 24	21 Sep 25	-
10	Hot Air Oven		Memmert	UF55 / B216.1666	National Food Institute, Ministry of Industry, Thailand	2500116-001-01	8 Oct 24	7 Oct 25	-
11	BOD Incubator	บีโอดี (Biochemical Oxygen Demand)	Arco	UC4-1320	Technology Promotion Association (Thailand-Japan)	24TJM1114	14 Jul 24	14 Jul 25	-
12	Water Bath	เบคเทอเรียลกลุ่มจีเอ็มแอล	Memmert	WNE 14 / L407.0756	Technology Promotion Association (Thailand-Japan)	24TJM1016/1	5 Aug 24	4 Aug 25	-
13	Autoclave		ALP	CL-40L / 807298	National Food Institute, Ministry of Industry, Thailand	2403982-001-01	7 Aug 24	6 Aug 25	-

List of Instruments Certification for Water Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration	Remark
Water									
14	Distillation Unit (Kjeldahl Method)	Total Kjeldahl Nitrogen (TKN) Ammonia	FOSS	KT9 Distillater / 91905393	FOSS South East Asia, Thailand	12875	5 Jul 24	5 Jul 25	-
15	Incubator	โคลิฟอร์มแบคทีเรียทั้งหมด (Total Coliform Bacteria) ฟิคอลโคลิฟอร์มแบคทีเรีย แบคทีเรียชนิดโคลิโด	BINDER	BD 53/ 13-07343	National Food Institute, Ministry of Industry, Thailand	2403982-003-01	7 Aug 24	6 Aug 25	-

CERTIFICATE OF ANALYSIS

Grade of Product: EPA PROTOCOL STANDARD

Customer: AIR LIQUIDE
(THAILAND) LTD
Part Number: E04NI99E15A01D3
Cylinder Number: EB0159156
Laboratory: 124 - Durham (SAP) - NC
PGVP Number: B22023
Gas Code: CO,NO,NOX,SO2,BALN

Reference Number: 122-402880224-1
Cylinder Volume: 144.0 CF
Cylinder Pressure: 2015 PSIG
Valve Outlet: 660
Certification Date: Nov 06, 2023

Expiration Date: Nov 06, 2026

Certification performed in accordance with "EPA Traceability Protocol for Assay and Certification of Gaseous Calibration Standards (May 2012)" document EPA 600/R-12/531, using the assay procedures listed. Analytical Methodology does not require correction for analytical interference. This cylinder has a total analytical uncertainty as stated below with a confidence level of 95%. There are no significant impurities which affect the use of this calibration mixture. All concentrations are on a mole/mole basis unless otherwise noted. The results relate only to the items tested. The report shall not be reproduced except in full without approval of the laboratory. Do Not Use This Cylinder below 100 psig, i.e. 0.7 megapascals.

ANALYTICAL RESULTS				
Component	Requested Concentration	Actual Concentration	Protocol Method	Total Relative Uncertainty
NOX	45.00 PPM	46.78 PPM	G1	+/- 1.3% NIST Traceable
NITRIC OXIDE	45.00 PPM	46.77 PPM	G1	+/- 1.3% NIST Traceable
SULFUR DIOXIDE	45.00 PPM	42.89 PPM	G1	+/- 1.0% NIST Traceable
CARBON MONOXIDE	1000 PPM	965.9 PPM	G1	+/- 0.7% NIST Traceable
NITROGEN	Balance			

CALIBRATION STANDARDS				
Type	Lot ID	Cylinder No	Concentration	Uncertainty
NTRM	21060703	CC707954	48.41 PPM NITRIC OXIDE/NITROGEN	+/- 1.2%
PRM	PRM	D913660	15.01 PPM NITROGEN DIOXIDE/AIR	+/- 1.5%
GMIS	124206889	CC322689	4.573 PPM NITROGEN DIOXIDE/NITROGEN	+/- 1.6%
NTRM	16061005	CC473180	49.02 PPM SULFUR DIOXIDE/NITROGEN	+/- 0.8%
NTRM	14060154	CC436951	990.9 PPM CARBON MONOXIDE/NITROGEN	+/- 0.6%

ANALYTICAL EQUIPMENT		
Instrument/Make/Model	Analytical Principle	Last Multipoint Calibration
Nicolet iS50 AUP2010249 CO	FTIR	Oct 11, 2023
Nicolet iS50 AUP2010249 NO	FTIR	Oct 11, 2023
Nicolet iS50 AUP2010249 NO2	FTIR	Oct 11, 2023
Nicolet iS50 AUP2010249 SO2	FTIR	Oct 11, 2023

Triad Data Available Upon Request

NOTES:GROSS WEIGHT: 28.4 kg

NET WEIGHT: 4.7 kg



Signature on file

Approved for Release

Page 1 of 1

CERTIFICATE OF ANALYSIS

Grade of Product: EPA Protocol

Part Number: E03AI99E15A006C
Cylinder Number: CC143232
Laboratory: 124 - Plumsteadville - PA
PGVP Number: A12020
Gas Code: CH4,PPN,BALA

Reference Number: 160-401908379-1
Cylinder Volume: 144.0 CF
Cylinder Pressure: 2016 PSIG
Valve Outlet: 590
Certification Date: Oct 16, 2020

Expiration Date: Oct 16, 2023

Certification performed in accordance with "EPA Traceability Protocol for Assay and Certification of Gaseous Calibration Standards (May 2012)" document EPA 600/R-12/531, using the assay procedures listed. Analytical Methodology does not require correction for analytical interference. This cylinder has a total analytical uncertainty as stated below with a confidence level of 95%. There are no significant impurities which affect the use of this calibration mixture. All concentrations are on a mole/mole basis unless otherwise noted. Do Not Use This Cylinder below 100 psig, i.e. 0.7 megapascals.

ANALYTICAL RESULTS				
Component	Requested Concentration	Actual Concentration	Protocol Method	Total Relative Uncertainty
METHANE	4000 PPM	4019 PPM	G1	+/- 1.0% NIST Traceable
PROPANE	4000 PPM	4008 PPM	G1	+/- 0.7% NIST Traceable
AIR	Balance			

CALIBRATION STANDARDS				
Type	Lot ID	Cylinder No	Concentration	Uncertainty
NTRM	02010405	K010090	4978 PPM PROPANE/NITROGEN	+/- 0.6%
NTRM	170608	CC160290	0.997 % METHANE/NITROGEN	+/- 0.4%

ANALYTICAL EQUIPMENT		
Instrument/Make/Model	Analytical Principle	Last Multipoint Calibration
MKS FTIR - CH4 - 000928781	FTIR	Oct 14, 2020
Nicolet 6700 APW1100391 C3H8	FTIR	Sep 18, 2020

Triad Data Available Upon Request

NOTES:NET WEIGHTS: 4.865kg

GROSS WEIGHTS: 27.365kg

PO#: 5220003825



Approved for Release

เอกสารไม่ควบคุม
Page 1 of 1

United Analyst and Engineering Consultant Co., Ltd.

3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260

Tel. 0 2763 2828 Fax 0 2763 2800 www.uaec consultant.com E-mail: uae@uaec consultant.com

United Analyst and Engineering Consultant Co., Ltd.

3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260

Tel. 0 2763 2828 Fax 0 2763 2800 www.uaec consultant.com E-mail: uae@uaec consultant.com

MULTI-POINT GAS TEST REPORT

Test Date : Sep 20, 2024

Equipment : Gas Analyzer (NO₂) Model : 42i
Manufacturer : Thermo Scientific Serial Number : CM22387040


Standard Gas Concentration		Dilutor Detail	
Sulphur Dioxide (SO ₂)	42.89 PPM	Manufacturer :	Thermo Scientific
Nitric Oxide (NO)	46.77 PPM	Model :	146i
Methane (CH ₄)	- PPM	Serial Number :	1180540071
Carbon Monoxide (CO)	965.9		
Cylinder No. :	EB0159156		
Expiration Date :	Nov 6, 2026		

Multi-point gas test data

Level	Reference Value (ppb)	Analyzer Display (ppb)	Difference Error	Percent Error	[% Error]
Level 1	Zero	0.0	0.00	0.00	0.00
Level 2	20.00%	100.0	100.7	0.70	0.70
Level 3	40.00%	200.0	200.6	0.30	0.30
Level 4	60.00%	300.0	300.9	0.30	0.30
Level 5	80.00%	400.0	400.0	0.00	0.00

Remark : Measuring Range 500.0 ppb
:Acceptable Limit \pm 5%

Multi-Point Gas Test Chart

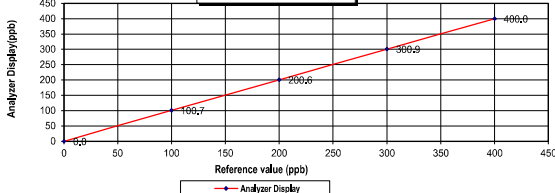

Calculate by
Gimhan.c
20/9/2567

Approve by
Pattom.k
20 Sep 2024

MULTI-POINT GAS TEST REPORT

Test Date : Sep 4, 2024

Equipment : Gas Analyzer (SO₂) Model : 43i
Manufacturer : Thermo SCIENTIFIC Serial Number : 1180540065


Standard Gas Concentration		Dilutor Detail	
Sulphur Dioxide (SO ₂)	42.89 PPM	Manufacturer :	Thermo SCIENTIFIC
Nitric Oxide (NO)	46.77 PPM	Model :	146i
Methane (CH ₄)	- PPM	Serial Number :	1180540071
Carbon Monoxide (CO)	965.9		
Cylinder No. :	EB01159156		
Expiration Date :	Nov 06, 2026		

Multi-point gas test data

Level	Reference Value (ppb)	Analyzer Display (ppb)	Difference Error	Percent Error	[% Error]
Level 1	Zero	0.0	0.00	0.00	0.00
Level 2	20.00%	100.0	100.7	0.70	0.70
Level 3	40.00%	200.0	200.4	0.20	0.20
Level 4	60.00%	300.0	300.8	0.27	0.27
Level 5	80.00%	400.0	400.0	0.00	0.00

Remark : Measuring Range 500.0 ppb
:Acceptable Limit \pm 5%

Multi-Point Gas Test Chart

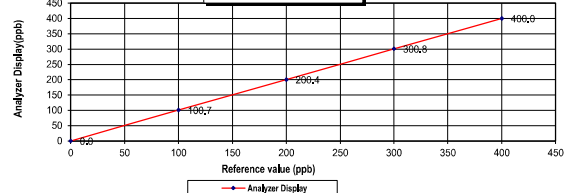

Calculate by
Gimhan.c
4/9/2567

Approve by
Pattom.k
4 Sep 2024

MULTI-POINT GAS TEST REPORT

Test Date : Sep 3, 2024

Equipment : Gas Analyzer (CO) Model : 48i
Manufacturer : Thermo Scientific Serial Number : CM08140004

Standard Gas Concentration

Sulphur Dioxide (SO₂) 42.89 PPM
Nitric Oxide (NO) 46.77 PPM
Methane (CH₄) - PPM
Carbon Monoxide (CO) 965.9 PPM
Cylinder No. : EB01159155
Expiration Date : Nov 06, 2026


Dilutor Detail

Manufacturer : Thermo Scientific
Model : 146i
Serial Number : 1180540071

Multi-point gas test data

	Reference Value (ppm)	Analyzer Display (ppm)	Difference Error	Percent Error	[% Error]
Level 1	Zero	0.0	0.0	0.0	0.0
Level 2	20.00%	10.0	10.5	0.5	4.8
Level 3	40.00%	20.0	20.9	0.9	4.3
Level 4	60.00%	30.0	30.7	0.7	2.3
Level 5	80.00%	40.0	40.0	0.0	0.0

Remark : Measuring Range 50.0 ppm
:Acceptable Limit $\pm 5\%$

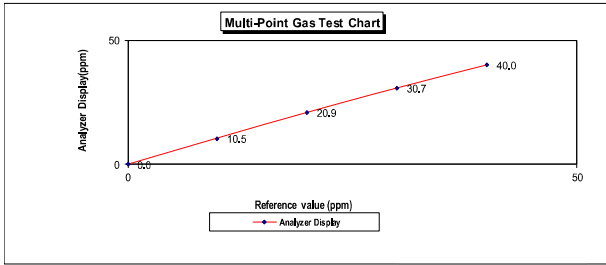

Calculate by
3...../.....9...../.....2567

Approve by
3...../.....9...../.....2567

Page 1 of 1

เอกสารไม่ควบคุม

MULTI-POINT GAS TEST REPORT

Test Date : Sep 2, 2024

Equipment : Hydrocarbon Analyzer Model : APHA-370
Manufacturer : HORIBA Serial Number : 93JN1MN9

Standard Gas Concentration

Sulphur Dioxide (SO₂) - PPM
Nitric Oxide (NO) - PPM
Methane (CH₄) 39.8 PPM
Carbon Monoxide (CO) - PPM
Cylinder No. : D824432
Expiration Date : Aug 4, 2028


Dilutor Detail

Manufacturer :
Model :
Serial Number :

Multi-point gas test data

	Reference Value (ppm)	Analyzer Display (ppm)	Difference Error	Percent Error	[% Error]
Level 1	Zero	0.00	0.37	0.37	0.37
Level 2	80.00%	40.00	39.56	-0.44	-1.11

Remark : Measuring Range 50.00 ppm
:Acceptable Limit $\pm 5\%$

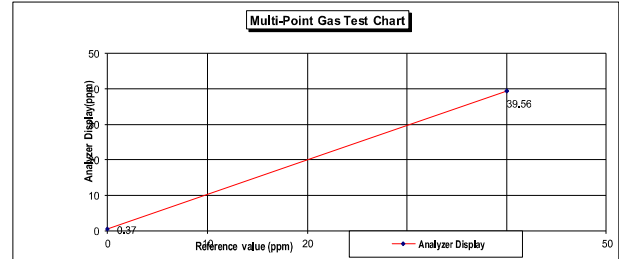


Calculate by
2...../.....9...../.....2567

Approve by
2...../.....9...../.....2567

Page 1 of 1

เอกสารไม่ควบคุม

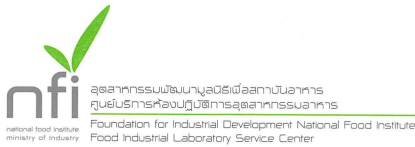

Calibration Certificate

Certificate No.: 2404042-001-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road,
Bangchack, Phrakhanong, Bangkok 10260

Page 1 of 5


Equipment: pH Meter
Manufacturer: HANNA Instruments
Model: HI2020-02
Serial No.: C0051107
ID No.: UAE WAO.005/2557
Order No.: 2404042
Operation No.: 2404042-001
Date of Receipt: 14 August 2024
Date of Calibration: 21 August 2024

Calibrated by Mr.Manas Somsak Specialist
Approved by (Mr.Pheraphat Tuanjit)
Manager, Division of Calibration Laboratory
Date of Issue: 21 August 2024
Responsible for the Technical Management Team

The uncertainties are for a confidence probability of approximately 95%.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65

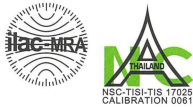

Calibration Report

Certificate No.: 2404042-001-01
Equipment: pH Meter Resolution: 0.001 pH ; 0.1 mV
Manufacturer: HANNA Instruments Model: HI2020-02
Serial No.: C0051107 Type: Bench top
ID No.: UAE.WAO.005/2557

Date of Calibration: 21 August 2024 Page 2 of 5

Location: Chemical Calibration Laboratory, National Food Institute
Environment Condition: Ambient Temperature: (23.1 \pm 1.5) °C Relative Humidity: (55 \pm 3) %
Condition of Equipment: Good Condition
Condition of this Results of Calibration

- Calibration Method
W-CC-002 : In house method based on direct measurement by using standard voltage calibrator and certified reference material (CRM)
- Reference Standards / Certified Reference Material

Instruments	Serial / ID No.	Manufacturer	Certificate No.	Due Date
2.1 DC Voltage Calibrator	2709007	Fluke	24E1752	30 May 2025
2.2 Digital Thermometer	2709007	Fluke	CC 660570-01	30 October 2024
2.3 Thermo-Hygro Meter	NFI.BTH 019/23	testo	QR24-0492	4 March 2025
Certified Reference Material				
2.4 pH buffer 4.008 (Primary pH buffer Solution)	873608	CPAchem	PH216.L5	16 February 2025
2.5 pH buffer 6.865 (Primary pH buffer Solution)	873609	CPAchem	PH217.L5	16 February 2025
2.6 pH buffer 10.01 (Primary pH buffer Solution)	949189	CPAchem	PH220.L5	30 November 2024
2.7 pH buffer 7.00 (Standard pH buffer Solution)	C03109	HACH LANGE GmbH	S11M004	16 October 2025
- This certification is traceable to The International System of Unit (SI Unit)

3.1 Instruments Ng.2.1	through	NSC-TISI-TIS 17025 Laboratory Accreditation of Calibration No.0008
3.2 Instruments Ng.2.2	through	NSC-TISI-TIS 17025 Laboratory Accreditation of Calibration No.0061
3.3 Instruments Ng.2.3	through	NSC-TISI-TIS 17025 Laboratory Accreditation of Calibration No.0292
3.4 Certified Reference Material Ng.2.4 to 2.6	traceable to	Primary measurement method: Hamed cell using calibrated thermometer, barometer, and nanovoltmeter. The Standard Solution preparation and certified by CPAchem Ltd is accredited to ISO 17034 and ISO/IEC 17025
3.5 Certified Reference Material Ng.2.7	traceable to	PTB Certificate Nr. PTB-PHQA-563/30504/23 and Certificate Nr. PTB-PHQB-555/30620/22 (PTB: Physikalisch-Technische Bundesanstalt, Braunschweig, Germany)

- This certificate was certified only for the instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.

F-CS-012 Revision: 01 Date: 20-04-65



Calibration Report

Certificate No.: 2404042-001-01
Equipment: pH Meter
Resolution: 0.001 pH ; 0.1 mV
Manufacturer: HANNA Instruments
Model: HI2020-02
Serial No.: C0051107
Type: Bench top
ID No.: UAE.WAO.005/2557

Date of Calibration: 21 August 2024 Page 3 of 5

Calibration Results:

1. Calibration of pH Meter (Manual Temperature Compensation at 25 °C)

Nominal pH	DC Voltage Standard (mV)	Average Indicator Reading		Uncertainty (±mV)	Coverage Factor (k)
		mV	pH		
0	414.122	421.0	-0.001	0.063	2.00
2	295.815	302.7	1.999	0.063	2.00
4	177.463	184.3	4.000	0.063	2.00
6	59.160	66.0	6.000	0.063	2.00
7	0.001	6.8	7.001	0.063	2.00
8	-59.159	-52.3	8.000	0.063	2.00
10	-177.462	-170.6	10.000	0.063	2.00
12	-295.813	-289.0	12.002	0.063	2.00
14	-414.121	-407.3	14.002	0.063	2.00

2. Calibration of pH Meter with Electrode (Manual Temperature Compensation at 25 °C)

Equipment: pH Electrode
Type: Combined Electrode
Manufacturer: HANNA Instruments
Model: HI11310
Serial No.: 539960
ID No.: N/A

Performance of Electrode system (Three-Point Calibration at pH 4, 7 and 10)

Certified Value @25 °C (pH)	Average Indicator Reading		Relative Slope (%)	Uncertainty (± pH)	Coverage Factor (k)
	pH	mV			
4.008	4.011	173.4	-	0.0043	2.00
7.001	7.004	-1.6	98.8	0.0073	2.00
9.997	10.011	-175.9	98.0	0.0073	2.00
6.865	6.870	6.4	-	0.0049	2.00

F-CS-012 Revision: 01 Date: 20-04-65

เอกสารไม่ควบคุม

Calibration Report

Certificate No.: 2404042-001-01
Equipment: Digital Thermometer with RTD (pH Meter)
Resolution: 0.1 °C
Model: HI2020-02
Serial No.: C0051107
ID No.: UAE.WAO.005/2557
Manufacturer: HANNA Instruments

Date of Calibration: 21 August 2024 Page 5 of 5

Calibration point: 15.0, 20.0 and 25.0 °C

Calibration result:

- The probe was immersed in liquid bath or dry bath to a minimum depth of 120 mm.
- Description of probe, model: HI11310 S/N: 539960
Dimension of probe : Diameter 12 mm., Length 120 mm.,
Sheath material : Glass

UUC* Reading (°C)	Standard Temperature (°C)	Correction Value (°C)	Uncertainty ± (°C)
15.0	14.998	0.0	0.099
20.0	19.999	0.0	0.099
25.0	24.999	0.0	0.099

Note

- UUC* : Unit Under Calibration

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k= 2, providing a level of confidence of approximately 95 %.

----- End -----

F-CS-012 Revision: 01 Date: 20-04-65

เอกสารไม่ควบคุม

Calibration Report

Certificate No.: 2404042-001-01
Equipment: Digital Thermometer with RTD (pH Meter)
Resolution: 0.1 °C
Model: HI2020-02
Serial No.: C0051107
ID No.: UAE.WAO.005/2557
Manufacturer: HANNA Instruments

Date of Calibration: 21 August 2024 Page 4 of 5

Location: Chemical Calibration Laboratory, National Food Institute
Environment Condition: Ambient Temperature 23 °C ± 1 °C
Relative Humidity 55 % ± 3 %

Condition of this results of Calibration:

- Calibration Method :
 - In house method: W-TE-025 by comparison with standard thermometer.
 - The Calibration is determined by comparing with a known temperature from a standard resistance thermometer.
 - The temperature scale in use at this laboratory is the International Temperature scale of 1990 (ITS-90).

2. Reference Standard Instrument :

Instrument	Model	Serial No.	Certificate No.	Due Date	Through
HANDHELD THERMOMETER	1521	A85997	TE 670101-01	16-Dec-24	NATIONAL FOOD INSTITUTE
Platinum Resistance Thermometer (PRT)	385	509201			

Support Equipment : - Low Temperature Bath (AMETEK RTC-187), Model: RTC-187C, S/N. 670930-00018

- This certificate is traceable to International System of Units (SI Units).
- This certificate was certified only for the instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.

6. Condition of Calibrated item : Good

7. Result of Calibration : ☒ Without adjustment ☐ After adjustment

F-CS-012 Revision: 01 Date: 20-04-65

เอกสารไม่ควบคุม

Calibration Certificate ID
TH2084-013-092024-ACC-TH

Mettler-Toledo (Thailand) Ltd.
846/4 - 846/5 Lesalle Rd., Bangna Tai Sub-District
Bangna District, Bangkok 10260
+662 723 0382
MT-TH.ServiceSupport@mt.com

METTLER TOLEDO

NSC-TISI-TIS 17025
CALIBRATION 0082

Accuracy Calibration Certificate

Customer

Company: United Analyst and Engineering Consultant Co.,Ltd.
Address: 3 Soi Udom Suk 41, Sukhumvit Rd., Bang Chak
City: Phra Khanong
Zip / Postal: 10260
State / Province: Bangkok
Order Number: 033323333

Weighing Device

Manufacturer: Mettler Toledo
Model: XPE205
Serial No.: B748058497
Building: N/A
Floor: 4
Room: Calibration Laboratory
Instrument Type: Weighing Instrument
Asset Number: UAE.CAL.004/2561
Terminal Model: PEAT
Terminal Serial No.: B748058497
Terminal Asset No.: N/A

Range	Max. Capacity	Readability (d)
1	220 g	0.00001 g

Procedure

Calibration Guideline: EURAMET cg-18 v. 4.0 (11/2015)
Mettler Toledo Work Instruction: CPW002/20

This calibration certificate contains measurements for As Found and As Left calibrations.

The sensitivity/span of the weighing instrument was adjusted before As Found and As Left calibrations with a built-in weight.
In accordance with EURAMET cg-18 (11/2015), the test loads were selected to reflect the specific use of the weighing device or to accommodate specific calibration conditions.

	Temperature		Humidity	
	As Found	Start: 22.1 °C End: 22.0 °C	Start: 45.8 % End: 43.4 %	
As Left	Start: 21.7 °C End: 22.3 °C	Start: 42.8 % End: 45.8 %		

As Found Calibration Date: 20-Sep-2024
As Left Calibration Date: 20-Sep-2024
Issue Date: 21-Sep-2024

Calibrator: Thanapong Insawang
Approved Signatory: Surachai P.
Technical Manager / Head of Calibration Center

Software Version: 1.24.0.306
Report Version: 2.19.19
Form Number: F1030

© METTLER TOLEDO
This is an original document and may not be partially reproduced without the written permission of the issuing calibration laboratory.

Page 1 of 6

เอกสารไม่ควบคุม

Error of Indication

As Found

	Tare Load	Reference Value	Indication	Error of Indication	Expanded Uncertainty	k
1	N/A	0.00000 g	0.00000 g	0.00000 g	0.029 mg	2
2	N/A	10.00001 g	10.00002 g	0.00001 g	0.065 mg	2
3	N/A	29.99998 g	30.00002 g	0.00004 g	0.14 mg	2
4 *	N/A	50.00000 g	50.00010 g	0.00010 g	0.12 mg	2
5 *	50 g	50.00000 g	50.00008 g	0.00008 g	0.12 mg	2
6 *	100 g	50.00000 g	50.00007 g	0.00007 g	0.12 mg	2
7 *	150 g	50.00000 g	50.00012 g	0.00012 g	0.12 mg	2
8 *	N/A	69.99997 g	70.00008 g	0.00011 g	0.20 mg	2
9	N/A	89.99999 g	90.00013 g	0.00014 g	0.27 mg	2
10 *	N/A	100.00006 g	100.00056 g	0.00050 g	0.20 mg	2
11 *	N/A	120.00003 g	120.00054 g	0.00051 g	0.31 mg	2
12 *	N/A	150.00006 g	150.00060 g	0.00054 g	0.31 mg	2
13	N/A	170.00003 g	170.00061 g	0.00058 g	0.38 mg	2
14 *	N/A	200.00025 g	200.00039 g	0.00014 g	0.35 mg	2

As Left

	Tare Load	Reference Value	Indication	Error of Indication	Expanded Uncertainty	k
1	N/A	0.00000 g	0.00000 g	0.00000 g	0.017 mg	2
2	N/A	10.00001 g	10.00002 g	0.00001 g	0.061 mg	2
3	N/A	29.99998 g	29.99997 g	-0.00001 g	0.14 mg	2
4 *	N/A	50.00000 g	50.00002 g	0.00002 g	0.12 mg	2
5 *	50 g	50.00000 g	50.00004 g	0.00004 g	0.12 mg	2
6 *	100 g	50.00000 g	50.00005 g	0.00005 g	0.12 mg	2
7 *	150 g	50.00000 g	50.00006 g	0.00006 g	0.12 mg	2
8 *	N/A	69.99997 g	69.99999 g	0.00002 g	0.20 mg	2
9	N/A	89.99999 g	90.00002 g	0.00003 g	0.27 mg	2
10 *	N/A	100.00006 g	100.00006 g	0.00000 g	0.20 mg	2
11 *	N/A	120.00003 g	120.00000 g	-0.00003 g	0.31 mg	2
12 *	N/A	150.00006 g	150.00001 g	-0.00005 g	0.31 mg	2
13	N/A	170.00003 g	170.00002 g	-0.00001 g	0.38 mg	2
14 *	N/A	200.00025 g	200.00027 g	0.00002 g	0.35 mg	2

*The calculated uncertainty was replaced by the CMC (Calibration and Measurement Capabilities) value because the calculated uncertainty was smaller than the CMC value.


Measurement Results

Repeatability

Test Load: 200 g

	As Found	As Left
1	200.00039 g	200.00023 g
2	200.00042 g	200.00023 g
3	200.00039 g	200.00022 g
4	200.00043 g	200.00023 g
5	200.00040 g	200.00022 g
6	200.00042 g	200.00021 g
7	200.00040 g	200.00023 g
8	200.00040 g	200.00022 g
9	200.00041 g	200.00022 g
10	200.00042 g	200.00021 g

Standard Deviation	0.000014 g	0.000008 g
--------------------	------------	------------

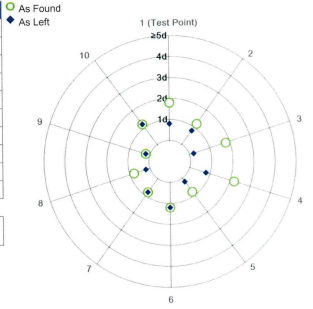


The "d" in the graph represents the readability of the range/interval in which the test was performed.
The results of this graph are based upon the absolute values of the differences from the mean value.

Eccentricity


Test Load: 100 g

Position	As Found	As Left
1	0.00000 g	0.00000 g
2	0.00000 g	0.00000 g
3	0.00000 g	0.00000 g
4	0.00000 g	0.00000 g
5	0.00000 g	0.00000 g

Maximum Deviation	0.00000 g	0.00000 g
-------------------	-----------	-----------

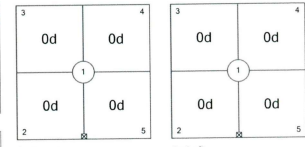

The "d" in the graph represents the readability of the range/interval in which the test was performed.

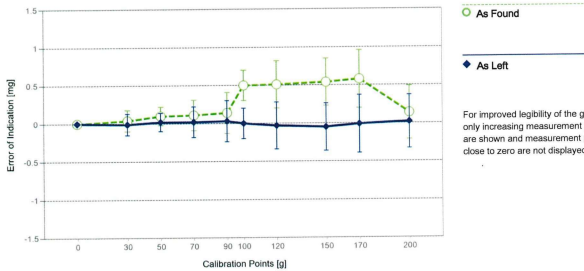

The expanded measurement uncertainty is reported as the standard measurement uncertainty multiplied by the coverage factor k such that the coverage probability corresponds to approximately 95 %.

The user is responsible for maintaining environmental conditions and the settings of the weighing instrument when it was calibrated.
The results of this calibration certificate relate only to the calibrated item.

Test Equipment

All weights used for metrological testing are traceable to national or international standards. The weights were calibrated and certified by an accredited calibration laboratory.

Weight Set 1: OIML E2

Weight Set No.:	WS01	Date of Issue:	31-May-2024
Certificate Number:	192393	Calibration Due Date:	30-Nov-2025

Weight Set 2: OIML E2

Weight Set No.:	WS01-1	Date of Issue:	18-Jun-2024
Certificate Number:	C428425034	Calibration Due Date:	17-Dec-2025

Weight Set 3: OIML E2

Weight Set No.:	WS01-2	Date of Issue:	12-Feb-2024
Certificate Number:	C408636201	Calibration Due Date:	11-Aug-2025

Thermo Hygrometer

Equipment No.:	IN278	Date of Issue:	11-Apr-2024
Certificate Number:	SG-H-00410/67	Calibration Due Date:	09-Apr-2025

Measurement Uncertainty of the Weighing Instrument in Use


Stated is the expanded uncertainty with $k=2$ in use. The formula shall be used for the estimation of the uncertainty under consideration of the errors of indication. The value R represents the net load indication in the unit of measure of the device.

Temperature coefficient for the evaluation of the measurement uncertainty in use: $1.0 \cdot 10^{-6} / K$
Temperature range on site for the evaluation of the measurement uncertainty in use: $3 K$

Range			As Found	As Left
d	Max			
1	0.00001 g	220 g	$U_1 = 0.029 \text{ mg} + 0.0322 \text{ mg/g} \cdot R$	$U_1 = 0.018 \text{ mg} + 0.00469 \text{ mg/g} \cdot R$

To optimize the stability of the linearization, besides of the zero load only increasing measurement points with a test load of 5% of the measurement range or larger are taken for the calculation of the linear equation.

Absolute and Relative Measurement Uncertainty in Use for Various Net Indications (Examples)				
Net Indication	As Found		As Left	
0.00220 g	0.029 mg	1.3%	0.018 mg	0.82%
0.02200 g	0.030 mg	0.14%	0.018 mg	0.082%
0.22000 g	0.036 mg	0.016%	0.019 mg	0.0087%
2.20000 g	0.100 mg	0.0045%	0.028 mg	0.0013%
220.00000 g	7.1 mg	0.0032%	1.0 mg	0.00048%

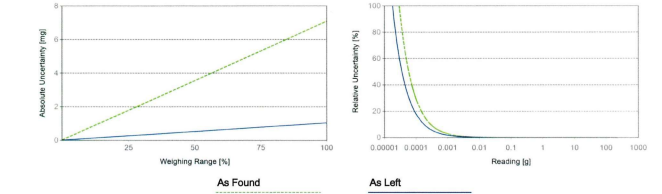

เอกสารไม่ควบคุม

Minimum Weight

As Found Minimum Weight Table

Minimum weights for different weighing tolerances and safety factors					
Tolerance	Safety Factor				
	1	2	3	5	10
0.1%	0.030104 g	0.062276 g	0.096738 g	0.173584 g	0.429425 g
0.2%	0.014806 g	0.030104 g	0.045918 g	0.079206 g	0.173584 g
0.5%	0.005865 g	0.011806 g	0.017825 g	0.030104 g	0.062276 g
1%	0.002923 g	0.005865 g	0.008826 g	0.014806 g	0.030104 g
2%	0.001459 g	0.002923 g	0.004392 g	0.007343 g	0.014806 g
5%	0.000583 g	0.001167 g	0.001752 g	0.002923 g	0.005865 g

✓ Pass: The determined minimum weight meets the requirement for the smallest net weight.

As Left Minimum Weight Table

Minimum weights for different weighing tolerances and safety factors					
Tolerance	Safety Factor				
	1	2	3	5	10
0.1%	0.017848 g	0.035696 g	0.054052 g	0.090953 g	0.186382 g
0.2%	0.008903 g	0.017848 g	0.026835 g	0.044937 g	0.090953 g
0.5%	0.003556 g	0.007111 g	0.010668 g	0.017848 g	0.035696 g
1%	0.001777 g	0.003556 g	0.005337 g	0.008903 g	0.017848 g
2%	0.000888 g	0.001777 g	0.002666 g	0.004444 g	0.008903 g
5%	0.000355 g	0.000711 g	0.001066 g	0.001777 g	0.003556 g

✓ Pass: The determined minimum weight meets the requirement for the smallest net weight.

At these net minimum weight values, the measurement uncertainty of the weighing device is equal to or less than 1/1 (no safety factor), 1/2, 1/3, 1/5, or 1/10 of the required tolerance. The values are calculated with $k = 2$ and based on the linear formula of the measurement uncertainty of the weighing device in use.

The safety factor for As Found is always 1. This implies no safety factor. As Found testing looks at the behavior of the instrument from the past until tested. For the past, it is necessary to know that the tolerance was met, but not the safety factor. The safety factor is a proactive measure to apply for future measurements.

Notes on minimum weight values in above table:

- If "N/A" is shown above, no appropriate value could be calculated.
- METTLER TOLEDO is not responsible for the definition of the process requirements.

เอกสารไม่ควบคุม

GWP®
Certificate



As Found



As Left



The weighing device meets the given process requirements.


The weighing device meets the given process requirements.

Tests Performed: ☒ As Found ☒ As Left

Process Requirements

Weighing Tolerance: 0.5% | Smallest Net Weight: 10.00000 g | Safety Factor: 2

Safe Weighing Range

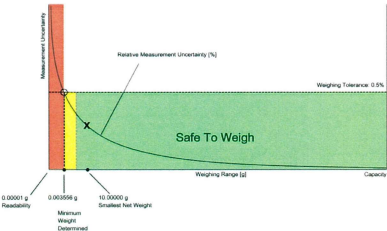

While the values in this graph reflect the actual calibration results, the measurement uncertainty curves are simply a visual representation. This graph reflects As Left testing, unless only As Found was performed.

เอกสารไม่ควบคุม

Measurement Results

Results Summary

	Repeatability	Eccentricity	Error of Indication
As Found	✓	✓	✓
As Left	✓	✓	✓

- ✓ = Passed
- ✗ = Failed
- ⚠ = Safety Factor not met

Repeatability

Test Load: 200 g

		As Found		As Left	
Tolerance	Control Limit	Std. Deviation	Result	Std. Deviation	Result
0.1%	0.005000 g	0.000014 g	✓	0.000008 g	✓
0.2%	0.010000 g		✓		✓
0.5%	0.025000 g		✓		✓
1%	0.050000 g		✓		✓
2%	0.100000 g		✓		✓
5%	0.250000 g		✓		✓

The weighing tolerance is met if the standard deviation is less than or equal to the corresponding control limit.

Eccentricity

Test Load: 100 g

		As Found		As Left	
Tolerance	Control Limit	Deviation	Result	Deviation	Result
0.1%	0.05000 g	0.00000 g	✓	0.00000 g	✓
0.2%	0.10000 g		✓		✓
0.5%	0.25000 g		✓		✓
1%	0.50000 g		✓		✓
2%	1.00000 g		✓		✓
5%	2.50000 g		✓		✓

The weighing tolerance is met if the deviation is less than or equal to the corresponding control limit.

เอกสารไม่ควบคุม

Error of Indication

As Found			Control limits for various weighing tolerances					
Tare	Reference Value	Error	0.1%	0.2%	0.5%	1%	2%	5%
N/A	0.00000 g	0.00000 g	N/A	N/A	N/A	N/A	N/A	N/A
N/A	29.99998 g	0.00004 g	0.01500 g	0.03000 g	0.07500 g	0.15000 g	0.30000 g	0.75000 g
N/A	50.00000 g	0.00010 g	0.02500 g	0.05000 g	0.12500 g	0.25000 g	0.50000 g	1.25000 g
50 g	50.00000 g	0.00008 g	0.02500 g	0.05000 g	0.12500 g	0.25000 g	0.50000 g	1.25000 g
100 g	50.00000 g	0.00007 g	0.02500 g	0.05000 g	0.12500 g	0.25000 g	0.50000 g	1.25000 g
150 g	50.00000 g	0.00012 g	0.02500 g	0.05000 g	0.12500 g	0.25000 g	0.50000 g	1.25000 g
N/A	69.99997 g	0.00011 g	0.03500 g	0.07000 g	0.17500 g	0.35000 g	0.70000 g	1.75000 g
N/A	89.99999 g	0.00014 g	0.04500 g	0.09000 g	0.22500 g	0.45000 g	0.90000 g	2.25000 g
N/A	100.00006 g	0.00050 g	0.05000 g	0.10000 g	0.25000 g	0.50000 g	1.00000 g	2.50000 g
N/A	120.00003 g	0.00051 g	0.06000 g	0.12000 g	0.30000 g	0.60000 g	1.20000 g	3.00000 g
N/A	150.00006 g	0.00054 g	0.07500 g	0.15000 g	0.37500 g	0.75000 g	1.50000 g	3.75000 g
N/A	170.00003 g	0.00058 g	0.08500 g	0.17000 g	0.42500 g	0.85000 g	1.70000 g	4.25000 g
N/A	200.00025 g	0.00014 g	0.10000 g	0.20000 g	0.50000 g	1.00000 g	2.00000 g	5.00000 g
Result			✓	✓	✓	✓	✓	✓

As Left			Control limits for various weighing tolerances					
Tare	Reference Value	Error	0.1%	0.2%	0.5%	1%	2%	5%
N/A	0.00000 g	0.00000 g	N/A	N/A	N/A	N/A	N/A	N/A
N/A	29.99998 g	-0.00001 g	0.01500 g	0.03000 g	0.07500 g	0.15000 g	0.30000 g	0.75000 g
N/A	50.00000 g	0.00002 g	0.02500 g	0.05000 g	0.12500 g	0.25000 g	0.50000 g	1.25000 g
50 g	50.00000 g	0.00004 g	0.02500 g	0.05000 g	0.12500 g	0.25000 g	0.50000 g	1.25000 g
100 g	50.00000 g	0.00005 g	0.02500 g	0.05000 g	0.12500 g	0.25000 g	0.50000 g	1.25000 g
150 g	50.00000 g	0.00006 g	0.02500 g	0.05000 g	0.12500 g	0.25000 g	0.50000 g	1.25000 g
N/A	69.99997 g	0.00002 g	0.03500 g	0.07000 g	0.17500 g	0.35000 g	0.70000 g	1.75000 g
N/A	89.99999 g	0.00003 g	0.04500 g	0.09000 g	0.22500 g	0.45000 g	0.90000 g	2.25000 g
N/A	100.00006 g	0.00000 g	0.05000 g	0.10000 g	0.25000 g	0.50000 g	1.00000 g	2.50000 g
N/A	120.00003 g	-0.00003 g	0.06000 g	0.12000 g	0.30000 g	0.60000 g	1.20000 g	3.00000 g
N/A	150.00006 g	-0.00005 g	0.07500 g	0.15000 g	0.37500 g	0.75000 g	1.50000 g	3.75000 g
N/A	170.00003 g	-0.00001 g	0.08500 g	0.17000 g	0.42500 g	0.85000 g	1.70000 g	4.25000 g
N/A	200.00025 g	0.00002 g	0.10000 g	0.20000 g	0.50000 g	1.00000 g	2.00000 g	5.00000 g
Result			✓	✓	✓	✓	✓	✓

The weighing tolerance is met if the error (of indication) for each test point is less than or equal to the corresponding control limit for that particular weighing tolerance. Results at or close to the zero point cannot be assessed.

Eccentricity

Test Load: 100 g		
Position	As Found	As Left
1	0.00000 g	0.00000 g
2	0.00000 g	0.00000 g
3	0.00000 g	0.00000 g
4	0.00000 g	0.00000 g
5	0.00000 g	0.00000 g
Maximum Deviation	0.00000 g	0.00000 g
Tolerance	0.000200 g ✓	0.000200 g ✓

The maximum deviation is determined as the absolute value of the largest deviation from the center.

Linearity - Differential Method

As Found

	Preload	Reference Value	Indication	Deviation
4	N/A	50.00000 g	50.00010 g	0.000008 g
5	50 g	50.00000 g	50.00008 g	-0.000005 g
6	100 g	50.00000 g	50.00007 g	-0.000028 g
7	150 g	50.00000 g	50.00012 g	0.000000 g
14*	N/A	200.00025 g	200.00039 g	N/A

Linearity Deviation	0.000028 g	Sensitivity Deviation	0.00014 g
Linearity Tolerance	0.0001 g ✓	Sensitivity Tolerance	0.0005 g ✓

The As Found Sensitivity Tolerance is only valid if the device has been adjusted before the test.

As Left

	Preload	Reference Value	Indication	Deviation
4	N/A	50.00000 g	50.00002 g	-0.000023 g
5	50 g	50.00000 g	50.00004 g	-0.000025 g
6	100 g	50.00000 g	50.00005 g	-0.000018 g
7	150 g	50.00000 g	50.00006 g	0.000000 g
14*	N/A	200.00025 g	200.00027 g	N/A

Linearity Deviation	0.000025 g	Sensitivity Deviation	0.00002 g
Linearity Tolerance	0.0001 g ✓	Sensitivity Tolerance	0.0005 g ✓

The values in column "Deviation" and the "Linearity Deviation" are zero point offset and sensitivity error compensated.

* This point was used to satisfy the sensitivity requirement.

Manufacturer Tolerance Assessment

Assessment done without considering measurement uncertainty.
The measurements from the attached calibration certificate were assessed against METTLER TOLEDO tolerances defined in SOP Test and Measurement Procedures for METTLER TOLEDO balances, Document: 10000018502.

	As Found	As Left
Overall	✓	✓
Repeatability	✓	✓
Eccentricity	✓	✓
Linearity	✓	✓
Sensitivity	✓	✓

Measurement Results

Repeatability

Test Load: 200 g		
	As Found	As Left
1	200.00039 g	200.00023 g
2	200.00042 g	200.00023 g
3	200.00039 g	200.00022 g
4	200.00043 g	200.00023 g
5	200.00040 g	200.00022 g
6	200.00042 g	200.00021 g
7	200.00040 g	200.00023 g
8	200.00040 g	200.00022 g
9	200.00041 g	200.00022 g
10	200.00042 g	200.00021 g
Standard Deviation	0.000014 g	0.000008 g
Tolerance	0.000030 g ✓	0.000030 g ✓

Preventive Maintenance for Laboratory Balances and
Halogen Moisture Analyzers

Customer Identification

Customer Name:	United Analyst and Engineering Consultant Co.,Ltd.		
Address:	3 Soi Udom Suk 41, Sukhumvit Rd., Bang Chak		
City:	Phra Khanong	State/Province:	Bangkok
Zip/Postal:	10260	Contact:	Suwit Chotnirk

Instrument Details

Manufacturer:	Mettler Toledo	Asset No.:	UAE.CAL.004/2561
Model:	XPE205	Alternate Asset No.:	N/A
Serial No.:	B748058497	Terminal Model:	PEAT
Building:	N/A	Terminal Serial No.:	B748058497
Floor:	4	Terminal Asset No.:	N/A
Room:	Calibration Laboratory	Customer Equip No.:	N/A
Instrument Type:	Single Range		

Procedure

Preventive Maintenance on this device is performed following the corresponding METTLER TOLEDO Work Instruction.

Service Date:	20-Sep-2024	Service Technician:	Thanapong Insawang
Next Service Date:	20-Sep-2025		
Date of Report Creation:	20-Sep-2024		

Instrument Settings

Software	Version
System SW Version	1.10 / 2.21

Adjustment Settings

Parameters	Setting
Adjustment Strategy	Internal
Trigger	Temp & Time
Schedule	Daily

Weighing Influencing Factors

Parameters	Level of Influence (High/Low/None)
Static Electricity	Low
Magnetism	None
Air Drafts	Low
Vibrations	Low
Heating Sources (e.g. Direct Sunlight)	None
Temperature Fluctuations	None

Cleaning

Task	Result
The instrument has been cleaned	Done

Visual Checks

Test	Result
Check housing	Pass
Check display	Pass
Check draft shield (if existing)	Pass
Check draft shield doors (if existing)	Pass
Check weighing pan	Pass
Check leveling	Pass
Check protective cover	Pass
Check AC power adapter	Pass
Check connectors and data cables	Pass
Check for correct fit and connection of accessories	Pass
Internal Inspection (if required)	N/A

Functional Tests

Test	Result
Perform startup of instrument	Pass
Perform touch screen test	Pass
Perform display test	Pass
Test draft shield motor movement (if existing)	Pass
Perform internal adjustment	Pass
Perform sensitivity test	Pass
Perform peripheral connectivity and function test	Pass
Perform SmartCal Test (only for Moisture Analyzers)	N/A
Perform Temperature Test (only for Moisture Analyzers)	N/A
Perform lid movement test (only for Moisture Analyzers)	N/A
Battery check	Pass
Check Date and Time	Pass

This document is issued to record completion of the work performed by METTLER TOLEDO on the subject device in accordance with agreed standards. It does not guarantee the continued performance of the subject device. Any measurements recorded are based on the subject device's performance at a given time as tested by METTLER TOLEDO and, except where explicitly stated otherwise, do not express an opinion as to the sufficiency of any customer designed procedures used to test the device. This document is not a warranty, either implied or express. METTLER TOLEDO expressly disclaims any liability arising from the use of the information in this document for any purpose other than as specified herein.

เอกสารไม่ควบคุม

Calibration Report

Certificate No.:	2500116-001-01
Equipment:	CHAMBER (Hot Air Oven)
Model:	UF55
Serial No.:	B216.1666
Resolution:	0.1 °C
ID No.:	UAE.WAO.027/2559
Manufacturer:	MEMMERT
Date of Calibration:	8 October 2024

Page 2 of 3

Location:	Laboratory, UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Environment Condition:	Ambient Temperature (30.3 ± 1) °C Relative Humidity (55 ± 1) % Line Voltage (230 ± 3) Volt

Condition of this results of Calibration:

- This instrument was calibrated by insert 9 standard thermometer into its chamber and calibration according to W-TE-014 Based on TLAS G-20-1/02-08 (E): Guidelines for Calibration and Checks of Temperature Controlled Enclosures.
 - The temperature scale used was based on ITS - 90.
 - All data show below were final values and the initial data may be obtained upon request.

2. Reference Standard Instrument :

Instrument	Model	Serial No./ID No.	Certificate No.	Due Date	Through
Digital Thermometer with sensor	34972A RTD	MYS7003188 CH#201-209/ RTD#201-209	TE 670486-01	8 June 2025	NATIONAL FOOD INSTITUTE

- This certificate is traceable to International System of Units (SI Units).
- This certificate was certified only for the instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.
- Condition of Calibrated item : Good

UUC Description :

Time of Record	1 Hour	9 Minute	At	104.0, 140.0 and 180.0 °C
Fresh air Damper	-	Open	Position	-
	X	Close	Fan	40%
	-	Not Available		

- Result of Calibration : ☒ Without adjustment ☐ After adjustment

เอกสารไม่ควบคุม

Calibration Certificate

Certificate No.:	2500116-001-01
Client name:	UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Address:	3 Soi Udomsuk 41, Sukhumvit Road, Bangchack, Prakhonong, Bangkok 10260

Page 1 of 3

Equipment: CHAMBER (Hot Air Oven)**Manufacturer:** MEMMERT**Model:** UF55**Serial No.:** B216.1666**ID No.:** UAE.WAO.027/2559**Order No.:** 2500116**Operation No.:** 2500116-001**Date of Receipt:** 8 October 2024**Date of Calibration:** 8 October 2024

Calibrated by	Mr.Yothin Charoensuk Scientist	Approved by	(Mr.Pheraphat Tuanjit) Manager, Division of Calibration Laboratory Responsible for the Technical Management Team
Date of Issue:	15 October 2024		

The uncertainties are for a confidence probability of approximately 95 %.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65

เอกสารไม่ควบคุม

Calibration Report

Certificate No.:	2500116-001-01
Equipment:	CHAMBER (Hot Air Oven)
Model:	UF55
Serial No.:	B216.1666
Resolution:	0.1 °C
ID No.:	UAE.WAO.027/2559
Manufacturer:	MEMMERT
Date of Calibration:	8 October 2024

Page 3 of 3

Calibration point: 104.0, 140.0 and 180.0 °C

Calibration Condition	Temperature (°C)	Relative Humidity (%)	Line Voltage (Volt)
MIN	29.3	54	227.0
MAX	31.2	56	232.0



Table 1 : Reporting of Temperature

Calibration point (°C)	Measured Temperature (°C) @ Sensor No. (Sensor No.9 is REF)									Uncertainty ± (°C)
	# 1	# 2	# 3	# 4	# 5	# 6	# 7	# 8	# 9	
104.0	103.89	103.66	103.88	103.89	104.40	103.98	103.70	104.10	104.15	0.53
140.0	139.85	139.53	139.87	139.88	140.67	140.00	139.60	140.25	140.23	0.73
180.0	179.63	179.22	179.71	179.76	181.03	180.06	179.41	180.87	180.39	0.90

Table 2 : Reporting of Characterization Result

UUC* Setting (°C)	UUC* Reading (°C)			Stability ± (°C)	Uniformity (°C)	Overall Variation (°C)
	MIN	MAX	Average			
104.0	104.0	104.0	104.0	0.15	0.49	0.88
140.0	140.0	140.0	140.0	0.13	0.71	1.2
180.0	180.0	180.0	180.0	0.13	1.2	1.9

Note The quoted uncertainty include " Stability " and " Loading effect (20% of Temp Uniformity) "

UUC* = Unit Under Calibration

Stability = One-half of the greatest maximum difference of measured temperatures at any one sensors, for at least half an hour after reaching steady state.

Uniformity = The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time.

Overall Variation = The difference of the maximum and minimum measured temperatures throughout observation time.

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k= 2, providing a level of confidence of approximately 95 %.

----- End -----

เอกสารไม่ควบคุม



Certificate of Calibration

Cert. No.: 24TM1113
Page : 1 of 3

Equipment : BOD Incubator
Manufacturer : ARCO
Model : UC4-1320
Serial No. : -
ID No. : UAE.WAO.002/2550
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260
Location : Lab Floor 2
Received Order : 11 July 2024
Calibration Date : 11 July 2024
Ambient Temperature : (26 ± 10) °C
Relative Humidity : (50 ± 30) %
Calibrated by : Tawatchai Pama
Approved by :
() Ponpan Palpim
(✓) Suwit Imjai
() Kunchit Promprat

Issue Date : 14 July 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2407-0243OC-1
Procedure Used :-

Cert. No.: 24TM1113
Page : 2 of 3

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD).
The temperature scale used was based on ITS-90.

Condition of this result of calibration

- Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Data Acquisition	MY49023932	23LM122	TPA	26 Jul 2024
- This certificate is valid only to the item calibrated on date and place of calibration.
- This certification is traceable to the International System of Unit.


Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Not Available

Environment during calibration		
	Beginning	Finished
Temp. (°C)	29	32
REL.Humid. (%)	78	65
AC Supply (Volt)	233	234

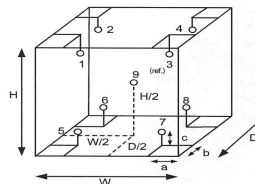

Probe Installation Details :

a = 10 cm
b = 10 cm
c = 10 cm

Dimension of Chamber :

D = 0.62 m
W = 1.2 m
H = 1.2 m
Capacity = 0.89 m³

Position :	Ref. Std. ID No.:
1	19-16RTD-01
2	19-16RTD-02
3	19-16RTD-03
4	19-16RTD-04
5	19-16RTD-05
6	19-16RTD-06
7	21-16RTD-07
8	19-16RTD-08
9 (ref.)	19-16RTD-09

เอกสารไม่ควบคุม



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2407-0243OC-1
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Not Available

Cert. No.: 24TM1113
Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor <i>k</i>
20.0	20.0	19.8	0.55	0.66	1.5	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	1	2	3	4	5	6	7	8	9 (ref.)	
20.0	20.210	20.331	20.162	19.645	20.287	20.070	19.838	19.781	19.954	0.79

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-000-

เอกสารไม่ควบคุม



Certificate of Calibration

Cert. No.: 24TM1016/1
Page : 1 of 3

This Certificate was issued to replace to the Certificate No. 24TM1016

Equipment : Water Bath

Manufacturer : Memmert

Model : WNB 14

Serial No. : L407.0756

ID No. : UAE.MIC.024/2550

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260

Location : Microbiology Laboratory

Received Order : 09 July 2024

Calibration Date : 09 - 10 July 2024

Ambient Temperature : (26 ± 10) °C

Relative Humidity : (50 ± 30) %

Calibrated by : Preecha Hlahib

Approved by :
Approved Signatory

() Ponpan Palpim
(✓) Suwit Imjai
() Kunchit Promprat

Issue Date : 5 August 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : Water Bath
Condition As-Received : Used Item
Reference : 2407-0153OC-6
Procedure Used :-

Cert. No.: 24TM1016/1
Page : 2 of 3

Calibration were conducted using in-house calibration procedure CP-OT04 Based on ASTM E715 according to direct measurement method with Data Acquisition which connected with Industrial Platinum Resistance Thermometer (IPRT).

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Data Acquisition	MY59003411	23LM208	TPA	27 Dec 2024

2. This certificate is valid only to the item calibrated on date and place of calibration.
3. This certification is traceable to the International System of Unit.


Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Heat transfer medium used : Water

	Environmental		AC Voltage Supply (Volt)
	(°C)	(%R.H.)	
Beginning of Calibration	26	67	221
Finished of Calibration	27	66	222

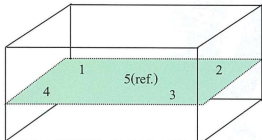

Front

Position :	Ref. Std. ID No.:
1	4804539-001
2	4804539-002
3	4804539-003
4	4804539-004
5(ref.)	4804539-005

เอกสารไม่ควบคุม

Calibration Certificate

Certificate No.: 2403982-001-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road,
Bangchack, Prakanhong, Bangkok 10260

Page 1 of 3

Equipment: Autoclave
Manufacturer: ALP
Model: CL-40L
Serial No.: 807298
ID No.: UAE.MIC.019/2560
Order No.: 2403982
Operation No.: 2403982-001
Date of Receipt: 7 August 2024
Date of Calibration: 7 August 2024

Calibrated by Mr.Manas Somsak Specialist
Approved by
(Miss Preeyaporn Jaengkarnkit)
Vice President, Department of Laboratory Services
Responsible for the Technical Management Team

Date of Issue: 14 August 2024

The uncertainties are for a confidence probability of approximately 95 %.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65



Equipment : Water Bath
Condition As-Received : Used Item
Reference : 2407-0153OC-6
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source

Cert. No.: 24TM1016/1
Page : 3 of 3

Calibration point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Average* Standard Reading (°C)					Uncertainty (± °C)
			Position					
			1	2	3	4	5(ref.)	
41.0	41.8	41.8	41.004	40.968	41.027	40.988	41.000	0.15
44.5	45.2	45.2	44.498	44.461	44.468	44.449	44.496	0.15
45.0	45.7	45.7	44.969	44.927	44.931	44.914	44.960	0.15

Calibration point (°C)	Uniformity (°C)	Stability (± °C)	Coverage Factor k
41.0	0.085	0.041	2
44.5	0.089	0.038	2
45.0	0.095	0.036	2

Average* : The average of 30 values in each position.

Uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Stability : One-half of the greatest maximum difference of measured temperature at any one probe.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity.

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k , providing a level of confidence of approximately 95 %.

-o-o-

เอกสารไม่ควบคุม

Calibration Report

Certificate No.: 2403982-001-01
Equipment: Autoclave
Model: CL-40L Serial No.: 807298
Resolution: 1 °C ID No.: UAE.MIC.019/2560
Manufacturer: ALP


Date of Calibration: 7 August 2024

Page 3 of 3

Calibration point: 121 °C

Calibration result:

Calibration Condition	Temperature (°C)	Relative Humidity (%)	Line Voltage (Volt)
Min	28.0	55	224
Max	30.0	65	226

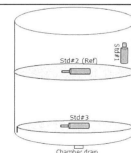

Standard at Position
Step1 = Attached to the load temperature probe, within 20 mm.
Step2 = In the upper half of the chamber.
Step3 = In the chamber drain, within 100 mm.

Table1 : Reporting of Temperature

Calibration Point (°C)	Measured Temperature (°C) @ Sensor No. (Sensor No.2 is REF)			Uncertainty (± °C)
	Std.# 1	Std.# 2 (Ref)	Std.# 3	
121	122.43	122.44	122.44	0.65

Table 2 : Reporting of Characterization Result

UUC* Setting (°C)	UUC* Reading			Stability (± °C)	Uniformity (°C)	Overall Variation (°C)
	Min (°C)	Max (°C)	Average (°C)			
122	122	122	122	0.11	0.065	0.031

Note

The quoted uncertainty include " Stability " and " Loading effect (20% of Uniformity)"

UUC* = Unit Under Calibration

Stability = One-half of the greatest maximum difference of measured temperatures at any one sensors, for at least half an hour after reaching steady state.

Uniformity = The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time.

Overall Variation = The Difference of the maximum and minimum measured temperatures throughout observation time.

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k = 2, providing a level of confidence of approximately 95 %.

----- End -----

F-CS-012 Revision: 01 Date: 20-04-65

Calibration Report

Certificate No.: 2403982-003-01
Equipment: CHAMBER (Incubator)
 Model: BD 53 Serial No.: 13-07343
 Resolution: 0.1 °C ID No.: UAE.MIC.005/2558
 Manufacturer: BINDER

Date of Calibration: 7 August 2024

Page 3 of 3

Calibration point: 44.0 °C

Calibration result:

Calibration Condition	Temperature (°C)	Relative Humidity (%)	Line Voltage (Volt)
MIN	24.8	50	227.0
MAX	24.9	60	232.0

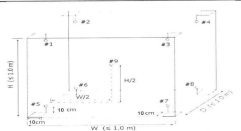


Table1 : Reporting of Temperature

Calibration point (°C)	Measured Temperature (°C) @ Sensor No. (Sensor No.9 is REF)									Uncertainty ± (°C)
	# 1	# 2	# 3	# 4	# 5	# 6	# 7	# 8	# 9	
44.0	44.24	44.33	44.30	44.46	43.42	43.51	43.40	43.53	43.97	0.30

Table 2 : Reporting of Characterization Result

UUC* Setting (°C)	UUC* Reading (°C)			Stability ± (°C)	Uniformity (°C)	Overall Variation (°C)
	MIN	MAX	Average			
44.0	44.0	44.0	44.0	0.083	0.57	1.2

Note: The quoted uncertainty include " Stability " and " Loading effect (20% of Temp Uniformity) "

UUC* = Unit Under Calibration

Stability = One-half of the greatest maximum difference of measured temperatures at any one sensors, for at least half an hour after reaching steady state.

Uniformity = The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time.

Overall Variation = The difference of the maximum and minimum measured temperatures throughout observation time.

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k= 2, providing a level of confidence of approximately 95 %.

----- End -----


 14 Aug. 2024

FCS-012 Revision: 01 Date: 20-04-65

2008 บางเขนพัฒนา 35 หมู่ 5 บางเขนพัฒนา แขวงบางเขน เขตบางเขน กรุงเทพมหานคร 10700

2008 Soi 35, Arun Amarin Road, Bang Yi Khan Subdistrict, Bang Phlat District, Bangkok 10700, Thailand

Tel: +66(0) 2422 6688 Fax: +66(0) 2422 8345

เอกสารไม่ลับ



nfi.or.th