

ภาคผนวก ง
เอกสารสอบเทียบเครื่องมือ

List of Instrument Certificates for Environmental Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*
1	Atomic Absorption Spectrometer	LEAD	Agilent Technologies	AA240FS / MY13160001	Agilent Technologies (Thailand) Co.,Ltd.	Preventive Maintenance Checklist	24/1/2024	23/1/2025
			Agilent Technologies	AA240FS / MY13160001	Agilent Technologies (Thailand) Co.,Ltd.	Preventive Maintenance Checklist	30/1/2025	29/1/2026
2	Analytical Balance	FAT OIL AND GREASE	Mettler Toledo	AB204-S/FACT / 1129361010	Technology Promotion Association (Thailand-Japan)	24MM292	11/5/2024	10/5/2025
			Mettler Toledo	AB204-S/FACT / 1129361010	United Analyst and Engineering Consultant Co., Ltd.	250422 1 BL002 25	23/4/2025	22/4/2026
3	Analytical Balance	TOTAL DISSOLVED SOLIDS TOTAL SOLIDS	Mettler Toledo	XSR205DU / C210685394	National Food Institute,Ministry of Industry, Thailand	2402283-002-01	2/4/2024	1/4/2025
			Mettler Toledo	XSR205DU / C210685394	National Food Institute,Ministry of Industry, Thailand	2502226-002-01	20/3/2025	19/3/2026
4	Analytical Balance	MIXED LIQUOR VOLATILE SUSPENDED SOLIDS SUSPENDED SOLIDS TOTAL SUSPENDED SOLIDS	Mettler Toledo	XSR205DU / C009071872	National Food Institute,Ministry of Industry, Thailand	2402283-001-01	2/4/2024	1/4/2025
			Mettler Toledo	XSR205DU / C009071872	National Food Institute,Ministry of Industry, Thailand	2502226-001-01	20/3/2025	19/3/2026
5	BOD Incubator	BIOCHEMICAL OXYGEN DEMAND	ARCO	UC4-1320 / 1021	Technology Promotion Association (Thailand-Japan)	24TM1113	11/7/2024	16/7/2025
6	BOD Incubator	BIOCHEMICAL OXYGEN DEMAND	ARCO	UC4-1320 / 13URC4S013201	Technology Promotion Association (Thailand-Japan)	24TM303	10/1/2024	8/2/2025
			ARCO	UC4-1320 / 13URC4S013201	Technology Promotion Association (Thailand-Japan)	25TM205	8/2/2025	7/2/2026

List of Instrument Certificates for Environmental Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*
7	BOD Incubator	BIOCHEMICAL OXYGEN DEMAND	ARCO	UC4-1320 / 1021	Technology Promotion Association (Thailand-Japan)	24TM1114	11/7/2024	10/7/2025
8	BOD Incubator	BIOCHEMICAL OXYGEN DEMAND	ARCO	UR-1320 / -	Technology Promotion Association (Thailand-Japan)	24TM588	1/4/2024	31/3/2025
9	BOD Incubator	BIOCHEMICAL OXYGEN DEMAND	ARCO	UR-1320 / -	Technology Promotion Association (Thailand-Japan)	24TM587	1/4/2024	31/3/2025
			ARCO	UR-1320 / -	Technology Promotion Association (Thailand-Japan)	25TM577	19/3/2025	19/3/2026
10	Continuous Flow Analyzer(CFA)	CYANIDE	Skalar Analytical B.V., the Netherlands	San++5000-02 / 182688	DKSH (Thailand) Ltd.	Service Report/Test Report WO-00018067	20/2/2024	19/2/2025
11	Cold Vapor Atomic Fluorescence Spectrometer (Mercury Analyzer)	MERCURY	Analytik Jena AG	Mercur duo plus / K 170A0153	Analytik Jena FarEast Thailand Ltd.	PM-OQ	3/2/2025	2/2/2026
12	DO Meter	DO	Horiba	LAQUA-DO210 / HE9M0028	Technology Promotion Association (Thailand-Japan)	24TW20	24/1/2024	22/1/2025
			Horiba	LAQUA-DO210 / HE9M0028	Technology Promotion Association (Thailand-Japan)	25TW19	23/1/2025	21/1/2026
13	DO Meter	DO	Horiba	LAQUA-DO210 / HE9M0021	Technology Promotion Association (Thailand-Japan)	24TW128	19/6/2024	17/6/2025
14	DO Meter	BIOCHEMICAL OXYGEN DEMAND	YSI	5100 / 11B 101863	Technology Promotion Association (Thailand-Japan)	24TW39	21/2/2024	20/2/2025
			YSI	5100 / 11B 101863	Technology Promotion Association (Thailand-Japan)	25TW29	18/2/2025	16/2/2026
15	SCT Meter	SALINITY	Horiba	LAQUA-EC210 / HC9L0012	Technology Promotion Association (Thailand-Japan)	24CH321	12/3/2024	10/3/2025
			Horiba	LAQUA-EC210 / HC9L0012	Technology Promotion Association (Thailand-Japan)	25CH244	26/2/2025	25/2/2026

List of Instrument Certificates for Environmental Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*
16	Gas Chromatography	TPH (C17-C35) TPH (C9-C16)	Agilent	GC 7890A / CN11021007	Agilent Technologies (Thailand) Co.,Ltd.	Certificate of System Qualification GC-OQ	18/2/2025	17/2/2026
17	Gas Chromatography/Mass Spectrometer	BENZENE n-HEXANE TOTAL XYLENES	Bruker	450GC / BR1201M099	Thai Unique Co.,Ltd	SV2407/21898	16/7/2024	15/7/2025
18	Hot Air Oven	TOTAL DISSOLVED SOLIDS TOTAL SUSPENDED SOLIDS	Memmert	UF55 / B212.0411	Technology Promotion Association (Thailand-Japan)	24TM589	1/4/2024	31/3/2025
19	Cooled Incubator	TOTAL PLATE COUNT	Binder	KB400 / WTB20200000015535	Technology Promotion Association (Thailand-Japan)	24TM647	1/4/2024	31/3/2025
			Binder	KB400 / WTB20200000015535	National Food Institute, Ministry of Industry, Thailand	2502229-006-01	19/3/2025	18/3/2026
20	Incubator	TOTAL PLATE COUNT	Binder	KB400 / 20220000022479	Technology Promotion Association (THAILAND-JAPAN)	24TM938	9/7/2024	8/7/2025
21	Incubator	TOTAL PLATE COUNT	Binder	KB400 / 20220000000391	Technology Promotion Association (Thailand-Japan)	24TM884	7/6/2024	6/6/2025
22	Incubator	LEGIONELLA SP.	Memmert	IPP260 / V616.0066	National Food Institute, Ministry of Industry, Thailand	2502229-002-01	19/3/2025	18/3/2026
23	Inductively Coupled Plasma- Optical Emission Spectrometer(ICP-OES)	LEAD	Agilent Technologies, USA	5110 VDV(G8015AA) / MY8030001	Agilent Technologies (Thailand) Co.,Ltd.	Preventive Maintenance Checklist	4/11/2024	3/11/2025
24	pH Meter and pH Electrode	pH	Mettler Toledo	pH S20 SevenEasyTM / 1231155210	National Food Institute Ministry of Industry, Thailand	2401718-001-01	11/3/2024	10/3/2025
			Mettler Toledo	pH S20 SevenEasyTM / 1231155210	National Food Institute Ministry of Industry, Thailand	2501844-001-01	24/2/2025	23/2/2026
25	pH Meter	pH	YSI Environmental	pH 100A / JC02729	Technology Promotion Association	24CH1070	27/8/2024	25/7/2025
26	SCT Meter	SALINITY	YSI Environmental	Pro 30 / 18K100976	Technology Promotion Association (Thailand-Japan)	24CH1381	7/11/2024	6/11/2025

List of Instrument Certificates for Environmental Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*
27	UV-VIS Spectrophotometer	PHENOLS CHEMICAL OXYGEN DEMAND	Hitachi	U-2900 / 21E22-009	DQE Services Co.,Ltd.	SP25-001	3/1/2025	2/1/2026
28	UV-VIS Spectrophotometer	SULPHIDE	Hitachi	U-2900 / 21E22-009	DQE Services Co.,Ltd.	SP24-001	4/1/2024	3/1/2025
			Hitachi	U-2900 / 21E22-009	DQE Services Co.,Ltd.	SP25-001	3/1/2025	2/1/2026
29	UV/VIS Spectrophotometer	CYANIDE NITRATE NITROGEN PHENOLS	Hitachi	U-5100 / 23A4-008	DQE Services Co.,Ltd.	SP24-028	11/9/2024	9/9/2025
30	Total Organic Carbon(TOC) Analyzer	TOTAL ORGANIC CARBON	Teledyne Tekmar	Lotix 15-1600-000 / US18038002	Saengvith 2000 Co.,Ltd.	MN-001-01-68-MN	6/2/2025	5/2/2026

Due Date of Calibration* : Based on the annual calibration plan. At least 1 time per year.



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000-29 FAX. 0-2719-9484



Cert.No.: 24CH321

Page.: 1 of 3

Certificate of Calibration

Equipment :	Conductivity Meter
Manufacturer :	Horiba
Model :	LAQUA-EC210
Serial No. :	HC9L0012
ID No. :	UAE.EFM.008/2563(EFM.SCT.02/63)
Condition As-Received:	Used Item
Received Date :	12 March 2024
Calibration Date :	14 March 2024
Reference :	2403-0387WSC-1
Submitted by :	United Analyst and Engineering Consultant Co.,Ltd. 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260
Ambient Temperature :	(25 ± 2.5) °C
Relative Humidity :	(50 ± 15) %
Calibration Procedure:	In -house method : - CP-CH6 by direct measurement with certified reference material (CRM) - CP-CH8 by Comparison with temperature standard

Calibrated by : Warakorn Lerngagtrakul

Approved by :

Approved Signatory

- () Pornthippa Tameyakul
() Unnopphol Harachai
(✓) Saithip Meangmai

Issue Date : 15 March 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

A 0064531



Cert.No.: 24CH321

Page.: 2 of 3

Condition of this result of calibration

1. Reference Standard Instrument :-

<u>Instrument</u>	<u>Serial No.</u>	<u>ID No.</u>	<u>Certificate No.</u>	<u>Due date</u>
1) Thermometer	1963878	130RC095	23I1051	05 Sep 2024
2) Ref. Std.Thermometer	4982054	110RC044	23I908	26 July 2024

- This Certification is traceable to SI Through Technology Promotion Association (Thailand - Japan)

2. Certified Reference Materials :-

- Conductivity calibration solution, CPA chem Ltd., The measurement results are traceable to SI through CPA chem Ltd., ANSI-ASQ National Accreditation Board, Accredited No. AR-1835

<u>Conductivity Solution</u>	<u>Manufacturer</u>	<u>Lot No.</u>	<u>Exp. date</u>
1413.0 $\mu\text{S/cm}$	CPA Chem	936624	19 Oct 2024
12.880 mS/cm	CPA Chem	931956	30 Sep 2024

- Control Conductivity calibration solution temperature by Water bath (25 ± 0.1) $^{\circ}\text{C}$

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration results

Function : Conductivity Measurement

(*) After Adjustment at 1413.0 $\mu\text{S/cm}$

Conductivity Electrode Serial No.: 9B9F0278

Standard Conductivity Solution	Before Adjustment UUC* Reading	After Adjustment UUC* Reading	Uncertainty of Measurement (\pm)	Coverage factor k
1413.0 $\mu\text{S/cm}$	1340 $\mu\text{S/cm}$	1413 $\mu\text{S/cm}$	9.2 $\mu\text{S/cm}$	2.00
12.880 mS/cm	13.70 mS/cm	14.27 mS/cm	0.086 mS/cm	2.00

Remark : - UUC* = Unit Under Calibration

Saitthip



Cert.No.: 24CH321

Page.: 3 of 3

Calibration Results

Function : Temperature Measurement

This equipment was connected with Temperature Probe;

- Model : 9383
- Serial No. : 9B9F0278

Dimension of probe;

- Length : 110 mm
- Diameter : 10 mm
- Immersion Depth : 90 mm

Calibration Result : Without adjustment

Calibration Point (°C)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty of Measurement (± °C)	Coverage factor <i>k</i>
25.0	25.002	25.0	-0.002	0.13	2.00
30.0	30.001	30.0	-0.001	0.13	2.00
35.0	35.002	35.0	-0.003	0.13	2.00

Remark : - UUC* = Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

Saidhmp



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert.No.: 24CH1070

Page.: 1 of 3

Equipment : pH Meter
Manufacturer : YSI
Model : pH100A
Serial No. : JC02729
ID No. : UAE.EFM.195/2561(ENV.pH.04/61)
Condition As-Received: Used Item
Received Date : 27 August 2024
Calibration Date : 28 August 2024
Reference : 2408-0882WSC-1
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong, Bangkok 10260

Ambient Temperature : (25 ± 2.5) °C
Relative Humidity : (50 ± 15) %
Calibration Procedure : In - house method :
- CP-CH5 by direct measurement with DC voltage
standard and direct measurement with
certified reference material (CRM)
- CP-CH8 by comparison with temperature standard

Calibrated by : Warakorn Lernqattrakul


Approved by : _____
Approved Signatory

() Unnopphol Harachai
() Ponpan Paipim
(✓) Saithip Meangmai

Issue Date : 29 August 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

 <input checked="" type="radio"/> PASS <input type="radio"/> NOT PASS	
Remarks: Temp: 0.5 (std thermometer) mv: 30 (std voltage input) pH: 0.05 (std buffer solution)	
1507 (1507) Verify	1507 (1507) Approve



Cert.No.: 24CH1070

Page.: 2 of 3

Condition of this calibration result

1. Reference Standard Instrument

<u>Instrument</u>	<u>Serial No.</u>	<u>ID No.</u>	<u>Cert. No.</u>	<u>Due Date</u>
1)Document Process Calibrator	43160066	130RC092	24E1320	22 Apr 2025
2)Ref. Standard Thermometer	2188080	130RC044	23I1216	10 Oct 2024

- This Certification is traceable to SI Through Technology Promotion Association (Thailand - Japan)

2. Certified Reference Materials :The measurement results are traceable to SI through Hach Lenge GmbH Li
Deutsche Akkreditierungsstelle, Accredited No.D-RM-15184-01-00

: The measurement results are traceable to SI through CPA chem Ltd.,
ANSI-ASQ National Accreditation Board, Accredited No. AR-1835

<u>Buffer Solution</u>	<u>Manufacturer</u>	<u>Lot No.</u>	<u>Exp. date</u>
pH 4.006	Hach Lenge GmbH	C03146	23 Feb 2026
pH 6.999	Hach Lenge GmbH	C03145	28 Feb 2026
pH 9.997	CPA chem	970853	25 Apr 2025

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration Results

Function : mV Measurement

Performing standard curve by Document Process Calibrator at pH (4,7)(7,10)

Unit Under Calibration	Nominal Value	Standard Voltage Input	Actual Reading		Uncertainty of Measurement (±mV)	Coverage factor <i>k</i>
	pH	mV	mV	pH		
pH Meter S/N.: JC02729	4.00	177.48	177	4.01	0.58	2.00
	7.00	0.00	0	7.00	0.58	2.00
	7.00	0.00	0	7.00	0.58	2.00
	10.00	-177.48	-177	10.01	0.58	2.00



Cert.No.: 24CH1070

Page.: 3 of 3

Calibration Results

Function : pH Measurement

Performing three buffers standard curve by using buffer nominal pH (4,7)(7,10)

Unit Under Calibration	Standard pH Buffer Solution	Actual pH Reading	Actual mV Reading (mV)	Uncertainty of pH Measurement (\pm)	Coverage factor k
pH Electrode S/N.: 231018SIA605377	4.006	4.01	173	0.0090	2.05
	6.999	7.00	-1	0.0084	2.00
	6.999	7.00	-1	0.0085	2.00
	9.997	10.00	-176	0.0092	2.00

Function : Temperature Measurement

(*) Without adjustment

This equipment was connected with Temperature Probe;

- Model : -

- Serial No. : 231018SIA605377

Dimension of probe

- Length : 110 mm.

- Diameter : 12 mm.

- Immersion Depth : 100 mm.

Calibration Point (°C)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty of measurement (\pm °C)	Coverage factor k
20.0	20.002	20.1	0.098	0.13	2.00
25.0	25.003	25.1	0.097	0.13	2.00
45.0	45.002	45.0	-0.002	0.13	2.00

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k , providing a level of confidence of approximately 95 %.

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert.No.: 24CH1381

Page.: 1 of 3

Equipment : Conductivity Meter
Manufacturer : YSI
Model : Pro 30
Serial No. : 18K100976
ID No. : UAE.EFM.071/2562(ENV.SCT.01/62)
Condition As-Received: Used Item
Received Date : 05 November 2024
Calibration Date : 06 November 2024
Reference : 2411-0124WSC-2
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road, Bangchak,
Phrakhanong, Bangkok 10260
Ambient Temperature : $(25 \pm 2.5) ^\circ\text{C}$
Relative Humidity : $(50 \pm 15) \%$
Calibration Procedure: In -house method :
- CP-CH6 by direct measurement
with certified reference material (CRM)
- CP-CH8 by comparison with temperature standard
Calibrated by : Warakorn Lerngagtrakul
Approved by :

Approved Signatory

☐ Unnopphol Harachai
☐ Ponpan Paipim
☒ Saithip Meangmai

Issue Date : 11 November 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Cert.No.: 24CH1381

Page.: 2 of 3

Condition of this result of calibration

1. Reference Standard Instrument :-

<u>Instrument</u>	<u>Serial No.</u>	<u>ID No.</u>	<u>Certificate No.</u>	<u>Due date</u>
1) Thermometer	1963878	130RC095	24I995	09 Sep 2025
2) Ref. Std.Thermometer	4982054	110RC044	24I757	14 July 2025

- This Certification is traceable to SI Through Technology Promotion Association (Thailand - Japan)

2. Certified Reference Materials :-

- Conductivity calibration solution, CPA chem Ltd., The measurement results are traceable to SI through CPA chem Ltd., ANSI-ASQ National Accreditation Board, Accredited No. AR-1835

<u>Conductivity Solution</u>	<u>Manufacturer</u>	<u>Lot No.</u>	<u>Exp. date</u>
1412.9 $\mu\text{S/cm}$	CPA Chem	1005307	15 June 2025
12.881 mS/cm	CPA Chem	1005308	15 June 2025

- Control Conductivity calibration solution temperature by Water bath (25 ± 0.1) $^{\circ}\text{C}$

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration results

Function : Conductivity Measurement

(*) After Adjustment at 1412.9 $\mu\text{S/cm}$

Conductivity Electrode Serial No.: 18L00008

Standard Conductivity Solution	Before Adjustment UUC* Reading	After Adjustment UUC* Reading	Uncertainty of Measurement (\pm)	Coverage factor k
1412.9 $\mu\text{S/cm}$	1146 $\mu\text{S/cm}$	1413 $\mu\text{S/cm}$	9.2 $\mu\text{S/cm}$	2.00
12.881 mS/cm	10.17 mS/cm	12.48 mS/cm	0.086 mS/cm	2.00

Remark : - UUC* = Unit Under Calibration



Cert.No.: 24CH1381

Page.: 3 of 3

Calibration Results

Function : Temperature Measurement

This equipment was connected with Temperature Probe;

- Model : 9383
- Serial No. : 18L00008

Dimension of probe;

- Length : 104 mm
- Diameter : 16 mm
- Immersion Depth : 90 mm

Calibration Result : Without adjustment

Calibration Point (°C)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty of Measurement (± °C)	Coverage factor <i>k</i>
15.0	15.004	14.9	-0.104	0.13	2.00
30.0	30.002	29.9	-0.102	0.13	2.00
45.0	45.005	45.0	-0.005	0.13	2.00

Remark : - UUC* = Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert.No.: 24MM292

Page.: 1 of 3

Equipment : Electronic Balance

Manufacturer : Mettler Toledo

Model : AB204-S/FACT

Serial No. : 1129361010

ID No. : UAE.WAS.002/2552

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260

Location : Balance Room (108)

Received order : 11 May 2024

Calibration Date : 11 May 2024

Ambient Temperature : 15 °C to 40 °C

Relative Humidity : 30 % to 90 %

Calibrated by : Khit Ruttanaprapachai

Approved by : 
Approved Signatory

() Ponpan Paipim
() Suwit Imjai
(✓) Kunchit Promprat

Issue Date : 15 May 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : Electronic Balance
Condition As-Received : Used Item
Reference : 2405-0166OC-1
Procedure used :-

Cert.No.: 24MM292

Page: 2 of 3

Calibration were conducted using in-house calibration procedure CP-OB01 based on UKAS LAB 14 according to direct measurement method against standard weight.

Condition of this result of calibration

1. Reference standard instruments:-

<u>Instruments</u>	<u>Model</u>	<u>Serial No.</u>	<u>ID No.</u>	<u>Test report No.</u>	<u>Due date</u>
1) Standard Weight Set (E2)	15884	24053	70RC007	MM-0013-24	25 Jan 2026

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This result of calibration was made on requested at the point specified by customer.

4. This certificate is not certified for any commercial transaction.

5. This certification is traceable to the International System of Unit.

Result of calibration () Without Adjustment (*) After Adjustment by Internal Calibration

Range capacity : 0 g to 220 g **Resolution** 0.0001 g

Before Adjustment :

<u>Applied Weight</u>	<u>Balance Reading</u>	<u>Correction</u>	<u>Measurement Uncertainty</u>	<u>Coverage Factor</u>
(g)	(g)	(g)	(± mg)	(k)
100	100.0000	0.0000	0.19	2.03
200	200.0006	-0.0006	0.30	2

After Adjustment :

1. **Determination of the standard deviation of weighing machine** (n = 10)

<u>Applied Weight</u>	<u>Standard Deviation of Reading (g)</u>
(g)	
100	0.00007
200	0.00005



Equipment : Electronic Balance
Condition As-Received : Used Item
Reference : 2405-0166OC-1


Cert.No.: 24MM292

Page: 3 of 3

Result of calibration

2. Effect of off center loading

A mass of 100 g was placed to various position on the pan.
The weighing machine reading error obtained is given in the table

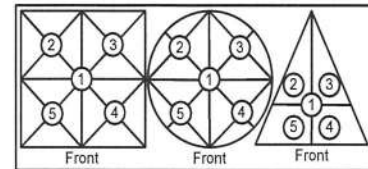

Maximum difference between
off-center and central loading

Position 1	Position 2	Position 3	Position 4	Position 5
(g)	(g)	(g)	(g)	(g)
-0.0004	-0.0004	-0.0003	-0.0003	-0.0004

(g)
0.0001

3. Departure from nominal value

Applied Weight	Balance Reading	Correction	Measurement Uncertainty	Coverage Factor
(g)	(g)	(g)	(\pm mg)	(k)
Unload	0.0000	0.0000	0.15	2.13
0.01	0.0100	0.0000	0.15	2.13
0.05	0.0500	0.0000	0.15	2.13
0.1	0.1000	0.0000	0.15	2.13
0.5	0.5000	0.0000	0.15	2.13
1	1.0000	0.0000	0.15	2.13
10	10.0000	0.0000	0.15	2.11
50	49.9999	+0.0001	0.17	2.06
100	99.9999	+0.0001	0.19	2.03
150	149.9998	+0.0002	0.29	2
200	199.9990	+0.0010	0.30	2

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k , providing a level of confidence of approximately 95 %.

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000-29 FAX. 0-2719-9484



Cert. No.: 24TM303

Page : 1 of 3

Certificate of Calibration

Equipment : BOD Incubator

Manufacturer : Arco

Model : UC4-1320

Serial No. : 13URC4S013201

ID No. : UAE.WAO.015/2561

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260

Location : Lab Floor 2

Received Order : 10 February 2024

Calibration Date : 10 February 2024

Ambient Temperature : (26 ± 10) °C

Relative Humidity : (50 ± 30) %

Calibrated by : Tawatchai Pama

Approved by :

Approved Signatory

- () Pornthippa Tameyakul
(☒) Unnopphol Harachai
() Suwit Imjai

Issue Date :

19 February 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2402-0234OC-1
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Not Available

Cert. No.: 24TM303

Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor <i>k</i>
20.0	20.1	19.9	0.37	0.72	1.4	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	
20.0	19.873	19.803	20.322	19.690	19.615	19.585	19.612	19.558	19.645	0.58

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2402-0234OC-1

Cert. No.: 24TM303

Page : 2 of 3

Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD).

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

<u>Instrument</u>	<u>Serial No.</u>	<u>Cert. No.</u>	<u>Traceable</u>	<u>Due Date</u>
1) Data Acquisition	MY59003411	23LM208	TPA	27 Dec 2024

2. This certificate is valid only to the item calibrated on date and place of calibration.


3. This certification is traceable to the International System of Unit.

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Not Available

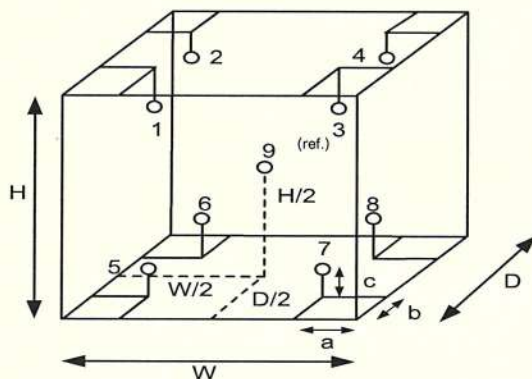

Environment during calibration		
	Beginning	Finished
Temp. (°C)	28	31
REL.Humid. (%)	70	65
AC Supply (Volt)	233	234

Position :	Ref. Std. ID No.:
1	20RTD-2/1
2	20RTD-2/2
3	20RTD-2/3
4	20RTD-2/4
5	20RTD-2/5
6	20RTD-2/6
7	20RTD-2/7
8	20RTD-2/8
9 (ref.)	20RTD-2/9

Probe Installation Details :

a = 10 cm
b = 10 cm
c = 10 cm

Dimension of Chamber :

D = 0.62 m
W = 1.2 m
H = 1.2 m
Capacity = 0.89 m³



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000-29 FAX. 0-2719-9484



Cert. No.: 24TM587
Page : 1 of 3

Certificate of Calibration

Equipment : BOD Incubator

Manufacturer : ARCO

Model : UR-1320

Serial No. : -

ID No. : UAE.WAO.018/2551

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260

Location : Lab Floor 2

Received Order : 01 April 2024
Calibration Date : 01 April 2024
Ambient Temperature : (26 ± 10) °C
Relative Humidity : (50 ± 30) %

Calibrated by : Krisda Malee

Approved by :

Approved Signatory

- () Ponpan Paipim
(☒) Suwit Imjai
() Kunchit Promprat

Issue Date : 5 April 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม

A 0065063



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2404-0004OC-1
Procedure Used :-

Cert. No.: 24TM587
Page : 2 of 3

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD).
The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

<u>Instrument</u>	<u>Serial No.</u>	<u>Cert. No.</u>	<u>Traceable</u>	<u>Due Date</u>
1) Data Acquisition	MY57013711	23LM115	TPA	11 Jul 2024

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

Remark : TPA : Technology Promotion Association (Thailand - Japan)


Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Not Available

Environment during calibration

	Beginning	Finished
Temp. (°C)	27	26
REL.Humid. (%)	48	49
AC Supply (Volt)	221	220

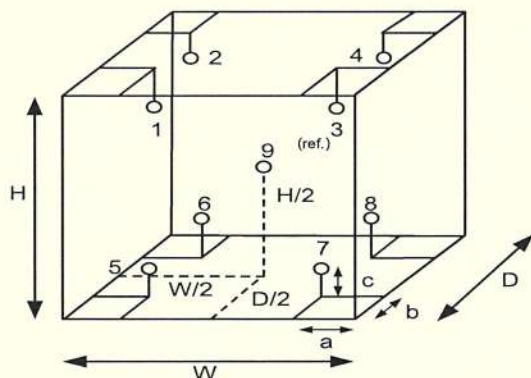

Position :	Ref. Std. ID No.:
1	18-18RTD-01
2	18-18RTD-02
3	18-18RTD-03
4	18-18RTD-04
5	18-18RTD-05
6	23-18RTD-06
7	18-18RTD-07
8	22-18RTD-08
9 (ref.)	18-18RTD-09

Probe Installation Details :

a = 10 cm
b = 10 cm
c = 10 cm

Dimension of Chamber :

D = 0.62 m
W = 1.2 m
H = 1.2 m
Capacity = 0.89 m³

Signature

เอกสารไม่ควบคุม
a 1209743



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2404-0004OC-1
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Not Available

Cert. No.: 24TM587

Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor <i>k</i>
20.0	20.0	20.0	0.45	0.55	1.3	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	
20.0	19.954	20.183	20.235	19.707	19.706	19.739	19.785	19.821	19.828	0.66

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่ควบคุม
a 1209742



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000-29 FAX. 0-2719-9484



Cert. No.: 24TM588

Page : 1 of 3

Certificate of Calibration

Equipment : BOD Incubator

Manufacturer : ARCO

Model : UR-1320

Serial No. : -

ID No. : UAE.WAO.006/2553

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260

Location : Lab Floor 2

Received Order : 01 April 2024

Calibration Date : 01 April 2024

Ambient Temperature : (26 ± 10) °C

Relative Humidity : (50 ± 30) %

Calibrated by : Krisda Malee

Approved by :

Approved Signatory

- () Ponpan Paipim
(☒) Suwit Imjai
() Kunchit Promprat

Issue Date :

5 April 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม

A 0065064



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2404-0004OC-2

Cert. No.: 24TM588

Page : 2 of 3

Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD).

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Data Acquisition	MY57013711	23LM115	TPA	11 Jul 2024

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

Remark : TPA : Technology Promotion Association (Thailand - Japan)


Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Not Available

Environment during calibration

	Beginning	Finished
Temp. (°C)	28	27
REL.Humid. (%)	45	47
AC Supply (Volt)	220	221

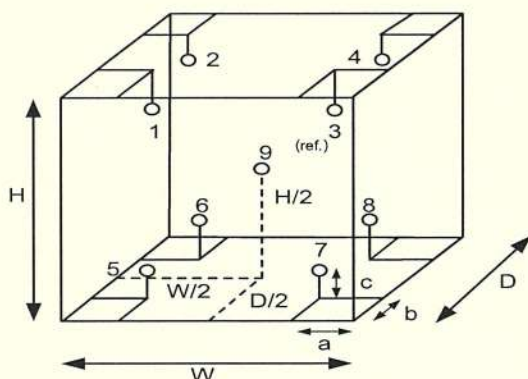

Position :	Ref. Std. ID No.:
1	22-18RTD-2/1
2	18RTD-2/2
3	18RTD-2/3
4	18RTD-2/4
5	18RTD-2/5
6	18RTD-2/6
7	18RTD-2/7
8	18RTD-2/8
9 (ref.)	18RTD-2/9

Probe Installation Details :

a = 10 cm
b = 10 cm
c = 10 cm

Dimension of Chamber :

D = 0.62 m
W = 1.2 m
H = 1.2 m
Capacity = 0.89 m³

เอกสารไม่ควบคุม

a 1209741



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2404-0004OC-2
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Not Available

Cert. No.: 24TM588

Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor <i>k</i>
20.0	20.0	19.9	0.47	0.69	1.4	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	
20.0	20.289	19.835	20.129	19.985	20.190	20.180	20.300	20.457	20.248	0.67

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่ควบคุม

a 1209740

กำหนดจุดห้ามใช้งาน

References Certificate Number. : 234TM588

Equipment : BOD Incubator

Model : UR-1320


Serial No. : -

ID No. : UAE.WAO.006/2553

Manufacturer : ARCO

Calibration Point : 20.0 °C

Unit Under Calibration Setting : 20.0 °C

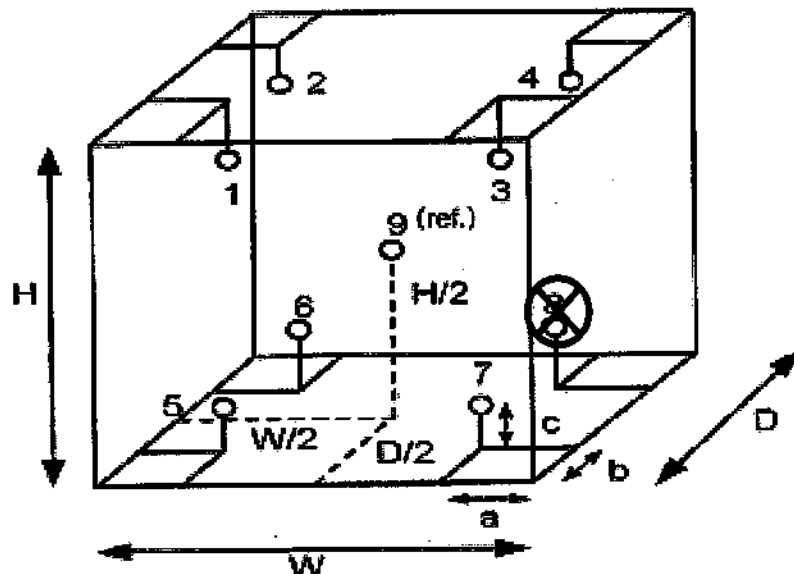

รูปภาพเครื่องมือ แสดงจุดที่ได้รับการสอบเทียบ และสัญลักษณ์ ⊗ แสดงจุดห้ามใช้งาน

กำหนดจุดห้ามใช้งานตำแหน่งที่....8.....

หมายเหตุ เก็บในแฟ้ม...../.....

\\uae.netapp\Lab-BK\INSTRUMENT (LI-2)\6.4\Certificate\ป้ายห้ามใช้งานเครื่องมือ\ป้ายห้ามใช้งานเครื่องมือ 2567\กำหนดจุดห้ามใช้งาน.doc

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000-29 FAX. 0-2719-9484



Cert. No.: 24TM589

Page : 1 of 3

Certificate of Calibration

Equipment : Hot Air Oven

Manufacturer : Memmert

Model : UF 55

Serial No. : B212.0411

ID No. : UAE.WAO.005/2556

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260

Location : Lab Floor 2

Received Order : 01 April 2024

Calibration Date : 01 - 02 April 2024

Ambient Temperature : (26 ± 10) °C

Relative Humidity : (50 ± 30) %

Calibrated by : Krisda Malee

Approved by :

Approved Signatory

- () Ponpan Paipim
(☒) Suwit Imjai
() Kunchit Promprat

Issue Date : 5 April 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม

A 0065065



Equipment : Hot Air Oven
Condition As-Received : Used Item
Reference : 2404-0004OC-3
Procedure Used :-

Cert. No.: 24TM589
Page : 2 of 3

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD) and Thermocouple Type T.

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-


<u>Instrument</u>	<u>Serial No.</u>	<u>Cert. No.</u>	<u>Traceable</u>	<u>Due Date</u>
1) Data Acquisition	MY57013711	23LM115	TPA	11 Jul 2024
2. This certificate is valid only to the item calibrated on date and place of calibration.				
3. This certification is traceable to the International System of Unit.				

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Close

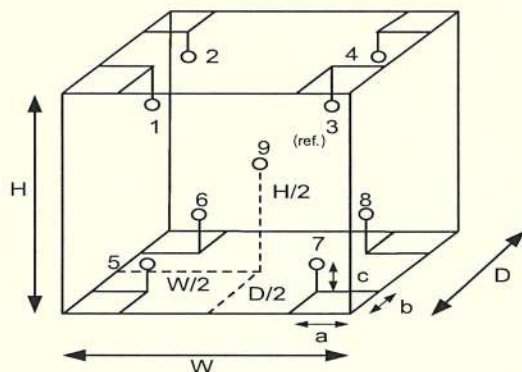

Environment during calibration		
	Beginning	Finished
Temp. (°C)	27	26
REL.Humid. (%)	47	48
AC Supply (Volt)	221	220

Ref. Std. ID No.: @ Calibration Point		
Position :	(120 to 180) °C	(104) °C
1	21-18TC-01	22-18RTD-2/1
2	21-18TC-02	18RTD-2/2
3	21-18TC-03	18RTD-2/3
4	21-18TC-04	18RTD-2/4
5	21-18TC-05	18RTD-2/5
6	21-18TC-06	18RTD-2/6
7	21-18TC-07	18RTD-2/7
8	21-18TC-08	18RTD-2/8
9 (ref.)	21-18TC-09	18RTD-2/9

Probe Installation Details : **Dimension of Chamber :**

a = 5.0 cm	D = 0.50 m
b = 5.0 cm	W = 0.80 m
c = 5.0 cm	H = 0.75 m
	Capacity = 0.30 m ³


เอกสารไม่ควบคุม
a 1209739



Equipment : Hot Air Oven
Condition As-Received : Used Item
Reference : 2404-0004OC-3
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Close

Cert. No.: 24TM589

Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor <i>k</i>
104.0	104.0	104.0	0.032	0.47	0.84	2
120.0	120.0	120.0	0.12	0.72	1.3	2
180.0	180.0	180.0	0.13	1.2	1.5	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (±°C)
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	
104.0	104.464	103.847	104.226	104.232	104.106	103.691	104.275	104.127	104.013	0.42
120.0	120.486	120.089	120.635	120.596	119.531	119.644	120.364	120.144	120.158	1.1
180.0	180.574	179.769	180.285	180.870	179.594	179.790	180.287	179.961	179.802	1.1

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่ควบคุม
a 1209738



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert. No.: 24TM647

Page : 1 of 3

Equipment : Incubator
Manufacturer : Binder
Model : KB 400 E6
Serial No. : 20200000015535
ID No. : UAE.MIC.018/2564
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260
Location : Microbiology Laboratory (302)
Received Order : 01 April 2024
Calibration Date : 01 April 2024
Ambient Temperature : (26 ± 10) °C
Relative Humidity : (50 ± 30) %

Calibrated by : Man Pattanapongpaiboon

Approved by :

Approved Signatory

() Ponpan Paipim
(✓) Suwit Imjai
() Kunchit Promprat

Issue Date : 7 April 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : Incubator
Condition As-Received : Used Item
Reference : 2404-0003OC-6

Cert. No.: 24TM647

Page : 2 of 3

Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD).

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

<u>Instrument</u>	<u>Serial No.</u>	<u>Cert. No.</u>	<u>Traceable</u>	<u>Due Date</u>
1) Data Acquisition	MY49023932	23LM122	TPA	26 Jul 2024

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.


Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Close

Environment during calibration		
	Beginning	Finished
Temp. (°C)	24	24
REL.Humid. (%)	54	57
AC Supply (Volt)	221	223

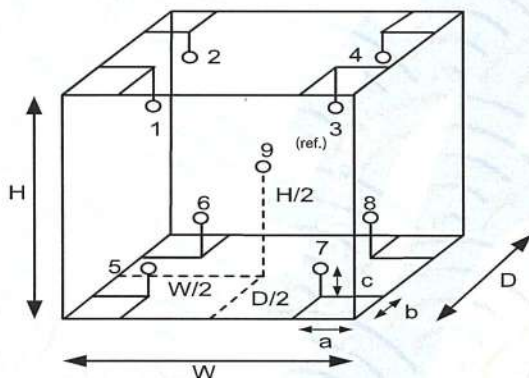

Position :	Ref. Std. ID No.:
1	20-16RTD-01
2	20-16RTD-02
3	20-16RTD-03
4	23-16RTD-04
5	22-16RTD-05
6	20-16RTD-06
7	20-16RTD-07
8	22-16RTD-08
9 (ref.)	22-16RTD-09

Probe Installation Details :

a = 10 cm
b = 10 cm
c = 10 cm

Dimension of Chamber :

D = 0.48 m
W = 0.65 m
H = 1.2 m
Capacity = 0.37 m³



Equipment : Incubator
Condition As-Received : Used Item
Reference : 2404-0003OC-6
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Close

Cert. No.: 24TM647

Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor <i>k</i>
35.0	35.0	35.0	0.035	0.19	0.22	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (±°C)
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	
35.0	35.000	35.022	34.841	34.851	35.027	35.011	35.023	35.028	35.007	0.30

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert. No.: 24TM884

Page : 1 of 3

Equipment : Incubator

Manufacturer : Binder

Model : KB 400

Serial No. : 20220000000391

ID No. : UAE.MIC.029/2565

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260

Location : Microbiology Laboratory


Received Order : 07 June 2024

Calibration Date : 07 June 2024

Ambient Temperature : (26 ± 10) °C

Relative Humidity : (50 ± 30) %

Calibrated by : Tawatchai Pama

Approved by : 
Approved Signatory

() Ponpan Paipim
() Suwit Imjai
(✓) Kunchit Promprat

Issue Date : 11 June 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : Incubator
Condition As-Received : Used Item
Reference : 2406-0190OC-2

Cert. No.: 24TM884

Page : 2 of 3

Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD).

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

<u>Instrument</u>	<u>Serial No.</u>	<u>Cert. No.</u>	<u>Traceable</u>	<u>Due Date</u>
1) Data Acquisition	MY49001451	24LM44	TPA	17 Mar 2025

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.


Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Close

Environment during calibration		
	Beginning	Finished
Temp. (°C)	21	19
REL.Humid. (%)	77	75
AC Supply (Volt)	228	229

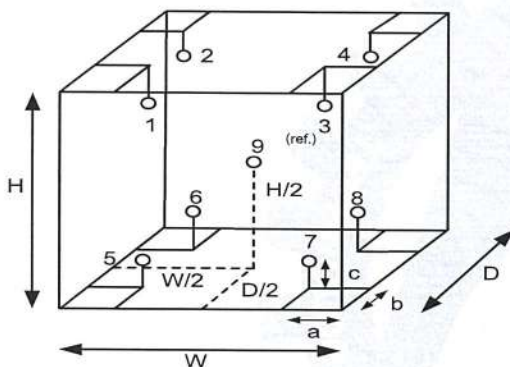

Position :	Ref. Std. ID No.:
1	19RTD-2/1
2	19RTD-2/2
3	19RTD-2/3
4	19RTD-2/4
5	19RTD-2/5
6	24-19RTD-2/6
7	19RTD-2/7
8	19RTD-2/8
9 (ref.)	19RTD-2/9

Probe Installation Details :

a = 10 cm
b = 10 cm
c = 10 cm

Dimension of Chamber :

D = 0.50 m
W = 0.65 m
H = 1.2 m
Capacity = 0.39 m³



Equipment : Incubator
Condition As-Received : Used Item
Reference : 2406-01900C-2
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Close

Cert. No.: 24TM884

Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor <i>k</i>
35.0	35.0	35.0	0.028	0.28	0.53	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	
35.0	35.317	35.184	35.142	35.064	35.098	35.093	34.894	34.826	35.056	0.30

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

กำหนดจุดห้ามใช้งาน

References Certificate Number. : 24TM884

Equipment : Incubator

Model : KB 400


Serial No. : 20220000000391

ID No. : UAE.MIC.029/2565

Manufacturer : Binder

Calibration Point : 35.0 °C

Unit Under Calibration Setting : 35.0 °C

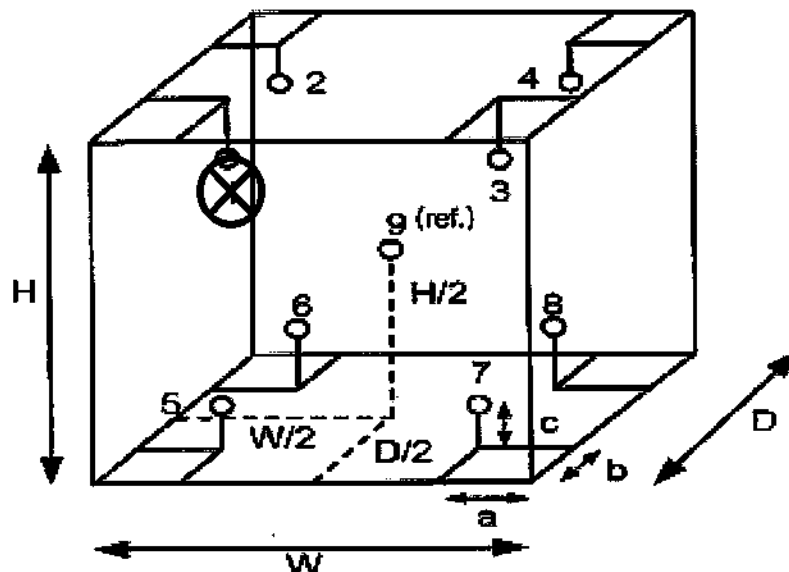

รูปภาพเครื่องมือ แสดงจุดที่ได้รับการสอบเทียบ และสัญลักษณ์ ⊗ แสดงจุดห้ามใช้งาน

กำหนดจุดห้ามใช้งานตำแหน่งที่....1.....

หมายเหตุ เก็บในแฟ้ม..... /

\\uae.netapp\Netapp_LAB\Lab-BK\INSTRUMENT (11-2)\6.4\Certificate\ป้ายห้ามใช้งานเครื่องมือ\ป้ายห้ามใช้งานเครื่องมือ 2567\กำหนดจุดห้ามใช้งาน

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert. No.: 24TM938

Page : 1 of 3

Equipment : Incubator
Manufacturer : Binder
Model : KB 400 E6
Serial No. : 20220000022479
ID No. : UAE.MIC.028/2566
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260
Location : Microbiology Laboratory
Received Order : 09 July 2024
Calibration Date : 09 July 2024
Ambient Temperature : (26 ± 10) °C
Relative Humidity : (50 ± 30) %

Calibrated by : Khit Ruttanaprapachai

Approved by :

Approved Signatory

- () Ponpan Paipim
(✓) Suwit Imjai
() Kunchit Promprat

Issue Date : 19 July 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : Incubator
Condition As-Received : Used Item
Reference : 2407-0153OC-4

Cert. No.: 24TM938

Page : 2 of 3

Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD).

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Data Acquisition	MY49001451	24LM44	TPA	17 Mar 2025

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.


Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Not Available

Environment during calibration		
	Beginning	Finished
Temp. (°C)	23	24
REL.Humid. (%)	52	54
AC Supply (Volt)	221	222

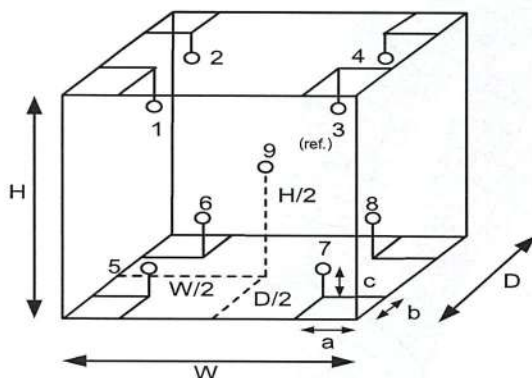

Position :	Ref. Std. ID No.:
1	19RTD-2/1
2	19RTD-2/2
3	19RTD-2/3
4	19RTD-2/4
5	19RTD-2/5
6	24-19RTD-2/6
7	19RTD-2/7
8	19RTD-2/8
9 (ref.)	19RTD-2/9

Probe Installation Details :

a = 10 cm
b = 10 cm
c = 10 cm

Dimension of Chamber :

D = 0.47 m
W = 0.65 m
H = 1.2 m
Capacity = 0.37 m³



Equipment : Incubator
Condition As-Received : Used Item
Reference : 2407-0153OC-4
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Not Available

Cert. No.: 24TM938

Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor <i>k</i>
35.0	35.0	35.0	0.030	0.31	0.33	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	
35.0	35.093	35.011	35.081	35.118	34.840	35.054	34.924	34.978	34.824	0.30

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k* , providing a level of confidence of approximately 95 %.

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert. No.: 24TM1113

Page : 1 of 3

Equipment : BOD Incubator

Manufacturer : ARCO

Model : UC4-1320

Serial No. : -

ID No. : UAE.WAO.002/2550

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260

Location : Lab Floor 2


Received Order : 11 July 2024

Calibration Date : 11 July 2024

Ambient Temperature : (26 ± 10) °C

Relative Humidity : (50 ± 30) %

Calibrated by : Tawatchai Pama

Approved by : 
Approved Signatory

() Ponpan Paipim
(✓) Suwit Imjai
() Kunchit Promprat

Issue Date : 14 July 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2407-0243OC-1

Cert. No.: 24TM1113

Page : 2 of 3

Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD).

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Data Acquisition	MY49023932	23LM122	TPA	26 Jul 2024

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.


Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Not Available

Environment during calibration		
	Beginning	Finished
Temp. (°C)	29	32
REL.Humid. (%)	78	65
AC Supply (Volt)	233	234

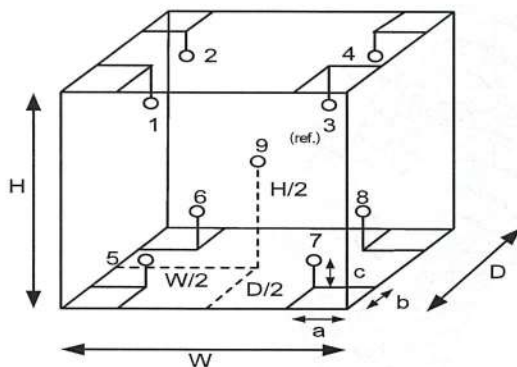

Position :	Ref. Std. ID No.:
1	19-16RTD-01
2	19-16RTD-02
3	19-16RTD-03
4	19-16RTD-04
5	19-16RTD-05
6	19-16RTD-06
7	21-16RTD-07
8	19-16RTD-08
9 (ref.)	19-16RTD-09

Probe Installation Details :

a = 10 cm
b = 10 cm
c = 10 cm

Dimension of Chamber :

D = 0.62 m
W = 1.2 m
H = 1.2 m
Capacity = 0.89 m³



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2407-0243OC-1
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Not Available

Cert. No.: 24TM1113

Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor <i>k</i>
20.0	20.0	19.8	0.55	0.66	1.5	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	
20.0	20.210	20.331	20.162	19.645	20.287	20.070	19.838	19.781	19.954	0.79

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert. No.: 24TM1114

Page : 1 of 3

Equipment : BOD Incubator

Manufacturer : ARCO

Model : UC4-1320

Serial No. : -

ID No. : UAE.WAO.018/2559

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260

Location : Lab Floor 2


Received Order : 11 July 2024

Calibration Date : 11 July 2024

Ambient Temperature : (26 ± 10) °C

Relative Humidity : (50 ± 30) %

Calibrated by : Tawatchai Pama

Approved by : 
Approved Signatory

() Ponpan Paipim
(✓) Suwit Imjai
() Kunchit Promprat

Issue Date : 14 July 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2407-0243OC-2

Cert. No.: 24TM1114

Page : 2 of 3

Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD).

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Data Acquisition	MY49023932	23LM122	TPA	26 Jul 2024

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.


Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Not Available

Environment during calibration		
	Beginning	Finished
Temp. (°C)	29	29
REL.Humid. (%)	78	72
AC Supply (Volt)	233	234

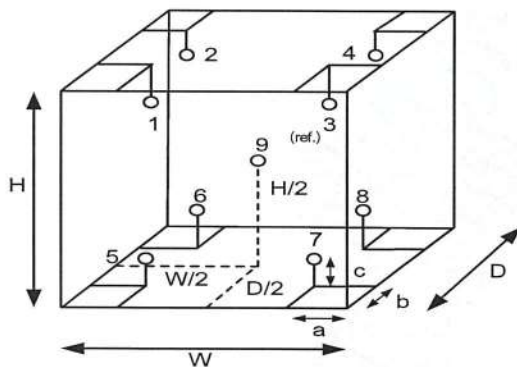

Position :	Ref. Std. ID No.:
1	20-16RTD-10
2	20-16RTD-02
3	20-16RTD-03
4	23-16RTD-04
5	22-16RTD-05
6	20-16RTD-06
7	20-16RTD-07
8	22-16RTD-08
9 (ref.)	22-16RTD-09

Probe Installation Details :

a = 10 cm
b = 10 cm
c = 10 cm

Dimension of Chamber :

D = 0.62 m
W = 1.2 m
H = 1.2 m
Capacity = 0.89 m³



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2407-0243OC-2
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Not Available

Cert. No.: 24TM1114

Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor <i>k</i>
20.0	20.0	19.9	0.29	0.81	1.2	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	
20.0	20.361	19.640	20.312	20.079	19.908	19.872	19.955	19.818	19.758	0.48

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3 : EQUIPMENT CALIBRATION AND TESTING SERVICES

534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250

TEL. 0-2717-3000 FAX. 0-2719-9484

Cert.No.: 24TW20

Page.: 1 of 2

Certificate of Testing

Equipment :	DO Meter
Manufacturer :	Horiba
Model :	LAQUA-DO210
Serial No. :	HE9M0028
ID No. :	UAE.EFM.013/2563 (EFM.DO.02/63)
Received Date :	23 January 2024
Test Date :	24 January 2024
Reference :	2401-0749WSC-1
Submitted by :	United Analyst and Engineering Consultant Co.,Ltd. 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260
Laboratory Condition :	Temperature (25 ± 5) °C Humidity (50 ± 20) %
Test Procedure :	In - house method : CP-CH9 by Comparison Technique with Azide Modification Method
Tested by :	Walalak Sirithean
Approved by :	 Approved Signatory
(<input checked="" type="checkbox"/>) Saithip Meangmai	
(<input type="checkbox"/>) Warakorn Lerngagtrakul	
(<input type="checkbox"/>) Ponpan Paipim	
Issue Date :	26 January 2024



Cert.No.: 24TW20

Page.: 2 of 2

Condition of this result of calibration

1. Reference Standard Instruments :

This certification is traceable to the International System of Unit through the reference standards laboratory of Industrial Calibration Center, Technology Promotion Association (Thailand-Japan).

<u>Instruments</u>	<u>Serial No.</u>	<u>ID No.</u>	<u>Certificate No.</u>	<u>Due Date</u>
1) Burette	-	130BU10	23CG1172	22 Mar 2025
2) Balance	1124013382	140RC006	23MM18	20 Feb 2024

2. Standard Material :-

<u>Material</u>	<u>Manufacturer</u>	<u>Lot.No.</u>	<u>Assay</u>
Sodium Thiosulfate pentahydrate	Merck	AM1763316	100.2%

Result : **Dissolved Oxygen Meter Adjustment With Air 100 %**

Dissolved Oxygen Probe No.: 9K9G0090

Titration Method (Azide Modification Method) (mg/L)	DO Meter Reading (mg/L)	Standard Deviation (mg/L)
8.20	8.20	0.0055

This report was certified only for the instrument we tested. It is allowable to use for study the system efficiency, The environmental impact control and present to organization it may concerned. Intend to use for advertising and referral purpose is prohibited. This report may not be reproduced other in full, without written approval of the laboratory

-o0o-

Saidhy

a 1199698



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000-29 FAX. 0-2719-9484



Cert. No.: 24LM9

Page.: 1 of 2

Certificate of Calibration

Equipment : DO Meter with Sensor

Manufacturer : Horiba

Model : LAQUA-DO210

Serial No. : HE9M0028

ID No. : UAE.EFM.013/2563(EFM.DO.02/63)

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260

Location : TPA Chemistry Calibration Laboratory

Received Order : 23 January 2024


Calibrated Date : 26 January 2024

Ambient Temperature : (26 ± 10) °C

Relative Humidity : (50 ± 30) %

AC Line Voltage : (220 ± 22) V

Calibrated by : Kunchit Promprat

Approved by : 
Approved Signatory

() Pornthippa Tameyakul
() Ponpan Paipim
(✓) Suwit Imjai

Issue Date : 31 January 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

A 0063154



Equipment : DO Meter with Sensor

Condition As-Received : Used Item

Reference : 2401-0749WSC-2

Cert. No.: 24LM9

Page.: 2 of 2

Procedure Used :-

Calibration were conducted using in-house calibration procedure CP-OT01 according to comparison with Industrial Platinum Resistance Thermometer (IPRT) into Temperature Bath.

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Digital Thermometer	3240076	231305	TPA	15 Mar 2024

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function : Temperature measurement.

This instrument was connected with temperature sensor, S/N.: 9K9G0090

Calibration Point (°C)	Immersion Depth (mm)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty (± °C)	Coverage Factor <i>k</i>
25.0	60	25.002	25.1	0.098	0.16	2.00
30.0	60	30.002	30.0	-0.002	0.16	2.00
35.0	60	35.001	35.0	-0.001	0.16	2.00

UUC* : Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

Yent



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3 : EQUIPMENT CALIBRATION AND TESTING SERVICES

534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250

TEL. 0-2717-3000 FAX. 0-2719-9484

Cert.No.: 24TW39

Page.: 1 of 2

Certificate of Testing

Equipment : DO Meter
Manufacturer : YSI
Model : 5100
Serial No. : 11B 101863
ID No. : UAE.WAO.004/2554
Received Date : 20 February 2024
Test Date : 21 February 2024
Reference : 2402-0629DSC-1
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road, Bangchak,
Phrakhanong, Bangkok 10260
Laboratory Condition : Temperature (25 ± 5) °C
Humidity (50 ± 20) %
Test Procedure : In - house method : CP-CH9
by Comparison Technique with Azide Modification Method

Tested by : Walalak Sirithean

Approved by :

Approved Signatory

- () Pornthippa Tameyakul
() Unnopphol Harachai
(☒) Saithip Meangmai

Issue Date : 22 February 2024

เอกสารไม่ควบคุม



Cert.No.: 24TW39

Page.: 2 of 2

Condition of this result of calibration

1. Reference Standard Instruments :

This certification is traceable to the International System of Unit through the reference standards laboratory of Industrial Calibration Center, Technology Promotion Association (Thailand-Japan).

<u>Instruments</u>	<u>Serial No.</u>	<u>ID No.</u>	<u>Certificate No.</u>	<u>Due Date</u>
1. Burette	-	130BU10	23CG1172	22 Mar 2025
2. Balance	14233821	110RC001	23MM405	16 July 2024

2. Standard Material :-

<u>Material</u>	<u>Manufacturer</u>	<u>Lot.No.</u>	<u>Assay</u>
Sodium Thiosulfate pentahydrate	Merck	AM1763316	100.2%

Result : Dissolved Oxygen Meter Adjustment With Air 100 %

Dissolved Oxygen Probe No.: 22B100125

Titration Method (Azide Modification Method) (mg/L)	DO Meter Reading (mg/L)	Standard Deviation (mg/L)
8.20	8.19	0.0055

This report was certified only for the instrument we tested. It is allowable to use for study
Intend to use for advertising and referral purpose is prohibited. This report may not be reproduced
other in full, without written approval of the laboratory

-o0o-

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3 : EQUIPMENT CALIBRATION AND TESTING SERVICES

534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250

TEL. 0-2717-3000 FAX. 0-2719-9484

Cert.No.: 24TW128

Page.: 1 of 2

Certificate of Testing

Equipment :	DO Meter
Manufacturer :	Horiba
Model :	LAQUA-DO210
Serial No. :	HE9M0021
ID No. :	UAE.EFM.014/2563 (EFM.DO.03/63)
Received Date :	18 June 2024
Test Date :	19 June 2024
Reference :	2406-0571WSC-1
Submitted by :	United Analyst and Engineering Consultant Co.,Ltd. 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260
Laboratory Condition :	Temperature (25 ± 5) °C Humidity (50 ± 20) %
Test Procedure :	In - house method : CP-CH9 by Comparison Technique with Azide Modification Method
Tested by :	Walalak Sirithean 
Approved by :	<hr/> <div>Approved Signatory</div>
() Unnopphol Harachai	
() Ponpan Paipim	
(✓) Saithip Meangmai	
Issue Date :	20 June 2024



Cert.No.: 24TW128

Page.: 2 of 2

Condition of this result of calibration

1. Reference Standard Instruments :

This certification is traceable to the International System of Unit through the reference standards laboratory of Industrial Calibration Center, Technology Promotion Association (Thailand-Japan).

<u>Instruments</u>	<u>Serial No.</u>	<u>ID No.</u>	<u>Certificate No.</u>	<u>Due Date</u>
1. Burette	-	130BU10	23CG1172	22 Mar 2025
2. Balance	14233821	110RC001	23MM405	16 July 2024

2. Standard Material :-

<u>Material</u>	<u>Manufacturer</u>	<u>Lot.No.</u>	<u>Assay</u>
Sodium Thiosulfate pentahydrate	Merck	AM1763316	100.2%

Result : **Dissolved Oxygen Meter Adjustment With Air 100 %**

Dissolved Oxygen Probe No.: 9K2H0067

Titration Method (Azide Modification Method) (mg/L)	DO Meter Reading (mg/L)	Standard Deviation (mg/L)
8.10	8.10	0.0071

This report was certified only for the instrument we tested. It is allowable to use for study
Intend to use for advertising and referral purpose is prohibited. This report may not be reproduced
other in full, without written approval of the laboratory

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert. No.: 24LM97

Page.: 1 of 2

Equipment : DO Meter with Sensor

Manufacturer : Horiba

Model : LAQUA-DO210

Serial No. : HE9M0021

ID No. : UAE.EFM.014/2563(EFM.DO.03/63)

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260

Location : TPA On Site Calibration Laboratory

Received Order : 18 June 2024


Calibrated Date : 20 June 2024

Ambient Temperature : (26 ± 10) °C

Relative Humidity : (50 ± 30) %

AC Line Voltage : (220 ± 22) V

Calibrated by : Warakorn Lerngagtrakul

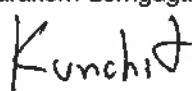Approved by : 
Approved Signatory

() Ponpan Paipim
() Suwit Imjai
(✓) Kunchit Promprat

Issue Date : 24 June 2024

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Equipment : DO Meter with Sensor
Condition As-Received : Used Item
Reference : 2406-0571WSC-2

Cert. No.: 24LM97
Page.: 2 of 2

Procedure Used :-

Calibration were conducted using in-house calibration procedure CP-OT01 according to comparison with Industrial Platinum Resistance Thermometer (IPRT) into Temperature Bath.

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

<u>Instrument</u>	<u>Serial No.</u>	<u>Cert. No.</u>	<u>Traceable</u>	<u>Due Date</u>
1) Digital Thermometer	2188080	2311216	TPA	11 Oct 2024

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function : Temperature measurement.

This instrument was connected with temperature sensor, S/N.: 9K2H0067

<u>Calibration Point</u> (°C)	<u>Immersion Depth</u> (mm)	<u>Standard Temperature</u> (°C)	<u>UUC* Reading</u> (°C)	<u>Error</u> (°C)	<u>Uncertainty</u> (± °C)	<u>Coverage Factor</u> <i>k</i>
25.0	80	25.002	25.0	-0.002	0.16	2.00
30.0	80	30.002	30.0	-0.002	0.16	2.00
35.0	80	35.003	35.0	-0.003	0.16	2.00

UUC* : Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert.No.: 25CH244

Page.: 1 of 3

Equipment : Conductivity Meter
Manufacturer : Horiba
Model : LAQUA-EC210
Serial No. : HC9L0012
ID No. : UAE.EFM.008/2563(EFM.SCT.02/63)
Condition As-Received: Used Item
Received Date : 25 February 2025
Calibration Date : 26 February 2025
Reference : 2502-0787WSC-1
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong, Bangkok 10260
Ambient Temperature : $(25 \pm 2.5) ^\circ\text{C}$
Relative Humidity : $(50 \pm 15) \%$
Calibration Procedure: In -house method :
- CP-CH6 by direct measurement
with certified reference material (CRM)
- CP-CH8 by comparison with temperature standard
Calibrated by : Warakorn Lerngagtrakul
Approved by :

Approved Signatory

() Chakrit Waewwanjua
() Ponpan Paipim
(✓) Saithip Meangmai

Issue Date : 27 February 2025

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Cert.No.: 25CH244

Page.: 2 of 3

Condition of this result of calibration

1. Reference Standard Instrument :-

<u>Instrument</u>	<u>Serial No.</u>	<u>ID No.</u>	<u>Certificate No.</u>	<u>Due date</u>
1) Thermometer	1963878	130RC095	24I995	09 Sep 2025
2) Ref. Std.Thermometer	4982054	110RC044	24I757	14 July 2025

- This Certification is traceable to SI Through Technology Promotion Association (Thailand - Japan)

2. Certified Reference Materials :-

- Conductivity calibration solution, CPA chem Ltd., The measurement results are traceable to SI through CPA chem Ltd., ANSI-ASQ National Accreditation Board, Accredited No. AR-1835

<u>Conductivity Solution</u>	<u>Manufacturer</u>	<u>Lot No.</u>	<u>Exp. date</u>
1412.9 $\mu\text{S/cm}$	CPA Chem	1005307	15 June 2025
12.881 mS/cm	CPA Chem	1005308	15 June 2025

- Control Conductivity calibration solution temperature by Water bath (25 ± 0.1) °C

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration results

Function : Conductivity Measurement

(*) After Adjustment at 1412.9 $\mu\text{S/cm}$

Conductivity Electrode Serial No.: 9B9F0278

Standard Conductivity Solution	Before Adjustment UUC* Reading	After Adjustment UUC* Reading	Uncertainty of Measurement (\pm)	Coverage factor <i>k</i>
1412.9 $\mu\text{S/cm}$	1390 $\mu\text{S/cm}$	1413 $\mu\text{S/cm}$	9.2 $\mu\text{S/cm}$	2.00
12.881 mS/cm	13.41 mS/cm	13.60 mS/cm	0.086 mS/cm	2.00

Remark : - UUC* = Unit Under Calibration



Cert.No.: 25CH244

Page.: 3 of 3

Calibration Results

Function : Temperature Measurement

This equipment was connected with Temperature Probe;

- Model : 9383
- Serial No. : 9B9F0278

Dimension of probe;

- Length : 110 mm
- Diameter : 16 mm
- Immersion Depth : 90 mm

Calibration Result : Without adjustment

Calibration Point (°C)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty of Measurement (± °C)	Coverage factor <i>k</i>
15.0	15.002	15.0	-0.002	0.13	2.00
30.0	30.002	30.0	-0.003	0.13	2.00
45.0	45.003	45.0	-0.003	0.13	2.00

Remark : - UUC* = Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-



Certificate of Calibration

Cert. No.: 25TM205

Page : 1 of 3

Equipment : BOD Incubator

Manufacturer : Arco

Model : UC4-1320

Serial No. : 13URC4S013201

ID No. : UAE.WAO.015/2561

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road
Bangchak, Phrakhanong
Bangkok 10260

Location : Lab. Floor 2

Received Order : 08 February 2025

Calibration Date : 08 February 2025

Ambient Temperature : (26 ± 10) °C

Relative Humidity : (50 ± 30) %

AC Line Voltage : (220 ± 22) V

Calibrated by :

Krisda Malee

Approved by :

Approved Signatory

() Chakrit Waewwanjua

() Suwit Imjai

(✓) Kunchit Promprat

Issue Date :

21 February 2025

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2502-0166OC-1

Cert. No.: 25TM205

Page : 2 of 3

Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD).

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

<u>Instrument</u>	<u>Serial No.</u>	<u>Cert. No.</u>	<u>Traceable</u>	<u>Due Date</u>
1) Data Acquisition	MY57013823	24LM71	TPA	12 May 2025

2. This certificate is valid only to the item calibrated on date and place of calibration.
3. This certification is traceable to the International System of Unit.


Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Not Available

Environment during calibration		
	Beginning	Finished
Temp. (°C)	26	25
REL.Humid. (%)	49	52
AC Supply (Volt)	221	220

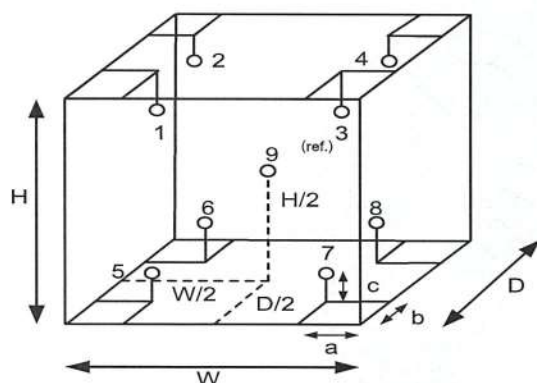

Position :	Ref. Std. ID No.:
1	21-17RTD-01
2	21-17RTD-02
3	17RTD-03
4	24-17RTD-04
5	17RTD-05
6	17RTD-06
7	17RTD-07
8	23-17RTD-08
9 (ref.)	23-17RTD-09

Probe Installation Details :

a = 10 cm
b = 10 cm
c = 10 cm

Dimension of Chamber :

D = 0.62 m
W = 1.2 m
H = 1.2 m
Capacity = 0.89 m³

เอกสารไม่ควบคุม



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2502-0166OC-1
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Not Available

Cert. No.: 25TM205

Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor <i>k</i>
20.0	20.0	19.9	0.36	0.56	0.99	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	(± °C)
20.0	19.841	19.714	20.110	19.862	19.747	19.710	19.676	19.789	19.695	0.54

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert. No.: 25TM577

Page : 1 of 3

Equipment : BOD Incubator

Manufacturer : ARCO

Model : UR-1320

Serial No. : -

ID No. : UAE.WAO.018/2551

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260

Location : Lab Floor 2

Received Order : 19 March 2025

Calibration Date : 19 March 2025

Ambient Temperature : $(26 \pm 10) ^\circ\text{C}$

Relative Humidity : $(50 \pm 30) \%$

AC Line Voltage : $(220 \pm 22) \text{ V}$

Calibrated by : Man Pattanapongpaiboon

Approved by :

Approved Signatory

- () Chakrit Waewwanjua
() Suwit Imjai
(✓) Kunchit Promprat

Issue Date : 27 March 2025

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2503-0437OC-1

Cert. No.: 25TM577

Page : 2 of 3

Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD).

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Data Acquisition	MY57013823	24LM71	TPA	12 May 2025

2. This certificate is valid only to the item calibrated on date and place of calibration.


3. This certification is traceable to the International System of Unit.

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Close

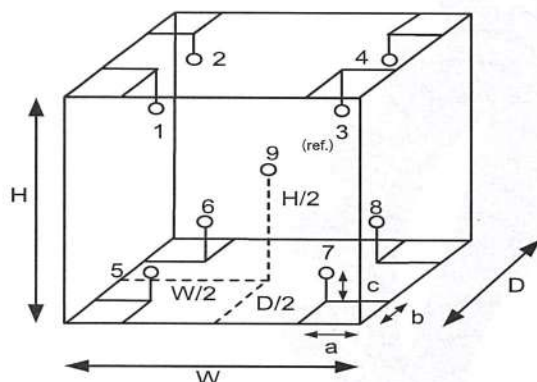

Environment during calibration		
	Beginning	Finished
Temp. (°C)	28	28
REL.Humid. (%)	56	55
AC Supply (Volt)	224	224

Position :	Ref. Std. ID No.:
1	21-17RTD-01
2	21-17RTD-02
3	17RTD-03
4	24-17RTD-04
5	17RTD-05
6	17RTD-06
7	17RTD-07
8	23-17RTD-08
9 (ref.)	23-17RTD-09

Probe Installation Details :

a = 10 cm
b = 10 cm
c = 10 cm

Dimension of Chamber :

D = 0.62 m
W = 1.2 m
H = 1.2 m
Capacity = 0.89 m³

เอกสารไม่ควบคุม



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2503-0437OC-1
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Close

Cert. No.: 25TM577

Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor <i>k</i>
20.0	20.0	20.0	0.24	0.54	0.99	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	
20.0	20.215	20.192	19.652	19.710	19.710	20.006	19.720	19.810	19.733	0.41

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3 : EQUIPMENT CALIBRATION AND TESTING SERVICES

534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250

TEL. 0-2717-3000 FAX. 0-2719-9484

Certificate of Testing

Cert.No.: 25TW19

Page.: 1 of 2

Equipment :	DO Meter
Manufacturer :	Horiba
Model :	LAQUA-DO210
Serial No. :	HE9M0028
ID No. :	UAE.EFM.013/2563 (EFM.DO.02/63)
Received Date :	22 January 2025
Test Date :	23 January 2025
Reference :	2501-0774WSC-1
Submitted by :	United Analyst and Engineering Consultant Co.,Ltd. 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260
Laboratory Condition :	Temperature (25 ± 5) °C Humidity (50 ± 20) %
Test Procedure :	In - house method : CP-CH9 by Comparison Technique with Azide Modification Method
Tested by :	Walalak Sirithean 
Approved by :	<hr/> <div>Approved Signatory</div>
() Pornthippa Tameyakul	
() Ponpan Paipim	
(✓) Saithip Meangmai	
Issue Date :	24 January 2025



Cert.No.: 25TW19

Page.: 2 of 2

Condition of this result of calibration

1. Reference Standard Instruments :

This certification is traceable to the International System of Unit through the reference standards laboratory of Industrial Calibration Center, Technology Promotion Association (Thailand-Japan).

<u>Instruments</u>	<u>Serial No.</u>	<u>ID No.</u>	<u>Certificate No.</u>	<u>Due Date</u>
1. Burette	-	130BU10	23CG1172	22 Mar 2025
2. Balance	14233821	110RC001	24MM131	04 July 2025

2. Standard Material :-

<u>Material</u>	<u>Manufacturer</u>	<u>Lot.No.</u>	<u>Assay</u>
Sodium Thiosulfate 5-Hydrate AR	KEMAUS	2203162447	99.6%

Result : Dissolved Oxygen Meter Adjustment With Air 100 %

Dissolved Oxygen Probe No.: 9K9G0090

Titration Method (Azide Modification Method) (mg/L)	DO Meter Reading (mg/L)	Standard Deviation (mg/L)
8.20	8.21	0.0045

This report was certified only for the instrument we tested. It is allowable to use for study
Intend to use for advertising and referral purpose is prohibited. This report may not be reproduced
other in full, without written approval of the laboratory

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert. No.: 25LM11

Page.: 1 of 2

Equipment : DO Meter with Sensor

Manufacturer : Horiba

Model : LAQUA-DO210

Serial No. : HE9M0028

ID No. : UAE.EFM.013/2563(EFM.DO.02/63)

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong, Bangkok 10260

Location : TPA On Site Calibration Laboratory

Received Order : 22 January 2025

Calibrated Date : 23 January 2025

Ambient Temperature : (26 ± 10) °C

Relative Humidity : (50 ± 30) %

AC Line Voltage : (220 ± 22) V

Calibrated by : Warakorn Lernagtrakul

Approved by :

Approved Signatory

- () Chakrit Waewwanjua
() Suwit Imjai
(✓) Kunchit Promprat

Issue Date : 29 January 2025

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Equipment : DO Meter with Sensor
Condition As-Received : Used Item
Reference : 2501-0774WSC-2
Procedure Used :-

Cert. No.: 25LM11
Page.: 2 of 2

Calibration were conducted using in-house calibration procedure CP-OT01 according to comparison with Industrial Platinum Resistance Thermometer (IPRT) into Temperature Bath.

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

<u>Instrument</u>	<u>Serial No.</u>	<u>Cert. No.</u>	<u>Traceable</u>	<u>Due Date</u>
1) Digital Thermometer	2188080	2411022	TPA	17 Sep 2025

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function : Temperature measurement.

This instrument was connected with temperature sensor, S/N.: 9K9G0090

<u>Calibration Point</u> (°C)	<u>Immersion Depth</u> (mm)	<u>Standard Temperature</u> (°C)	<u>UUC* Reading</u> (°C)	<u>Error</u> (°C)	<u>Uncertainty</u> (± °C)	<u>Coverage Factor</u> <i>k</i>
15.0	90	15.002	15.0	-0.002	0.16	2.00
30.0	90	30.002	30.0	-0.002	0.16	2.00
45.0	90	45.003	45.0	-0.003	0.16	2.00

UUC* : Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3 : EQUIPMENT CALIBRATION AND TESTING SERVICES


534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250

TEL. 0-2717-3000 FAX. 0-2719-9484

Certificate of Testing

Cert.No.: 25TW29

Page.: 1 of 2

Equipment : DO Meter
Manufacturer : YSI
Model : 5100
Serial No. : 11B 101863
ID No. : UAE.WAO.004/2554
Received Date : 14 February 2025
Test Date : 17 February 2025
Reference : 2502-0473DSC-1
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road, Bangchak,
Phrakhanong, Bangkok 10260
Laboratory Condition : Temperature (25 ± 5) °C
Humidity (50 ± 20) %
Test Procedure : In - house method : CP-CH9
by Comparison Technique with Azide Modification Method
Tested by : Walalak Sirithean

Approved by : _____
Approved Signatory
() Chakrit Waewwanjua
() Ponpan Paipim
(✓) Saithip Meangmai
Issue Date : 18 February 2025

เอกสารไม่ควบคุม



Cert.No.: 25TW29

Page.: 2 of 2

Condition of this result of calibration

1. Reference Standard Instruments :

This certification is traceable to the International System of Unit through the reference standards laboratory of Industrial Calibration Center, Technology Promotion Association (Thailand-Japan).

<u>Instruments</u>	<u>Serial No.</u>	<u>ID No.</u>	<u>Certificate No.</u>	<u>Due Date</u>
1. Burette	-	130BU10	23CG1172	22 Mar 2025
2. Balance	14233821	110RC001	24MM131	04 July 2025

2. Standard Material :-

<u>Material</u>	<u>Manufacturer</u>	<u>Lot.No.</u>	<u>Assay</u>
Sodium Thiosulfate 5-Hydrate AR	KEMAUS	2203162447	99.6%

Result : **Dissolved Oxygen Meter Adjustment With Air 100 %**
Dissolved Oxygen Probe No.: 24F100202

Titration Method (Azide Modification Method) (mg/L)	DO Meter Reading (mg/L)	Standard Deviation (mg/L)
8.22	8.22	0.0055

This report was certified only for the instrument we tested. It is allowable to use for study
Intend to use for advertising and referral purpose is prohibited. This report may not be reproduced
other in full, without written approval of the laboratory

-o0o-

เอกสารไม่ควบคุม

บันทึกผลการทวนสอบใบรับรองการสอบเทียบ (Verification of Certificate)

Certificate No. : 25TW29				Equipment : Do Meter			
Brand : YSI				Model : 5100			
Serial No. : 11B 101863				ID No. : UAE.WAO.004/2554			
Calibration results							
Titration Method	Standart Deviation	Do meter Reading	Error%	Correction%	Error Total Error	Judgement	(Total Error < Judgement)
(mg/L)	(mg/L)	(mg/L)	(mg/L)	(mg/L)	(mg/L)	(± mg/L)	(mg/L)
8.22	0.0055	8.22	0.0000	0.0000	0.0	0.02	pass
ผู้บันทึก.....อิสรา.บุญประกอบ.....ผู้ตรวจสอบ..... พ.รังษวิมล							
วันที่.....28/02/2025.....วันที่..... 28 ก.พ. 68							
หมายเหตุ :							

เก็บใบนี้เพื่อ.....

...../.....

เอกสารไม่ควบคุม

Certificate No.: 250422-1-BL002-25

Code No.: BL002-25

Page: 2 of 3

Equipment: Electronic Balance

Manufacturer: Mettler Toledo

Model: AB204-S/FACT

Readability: 0.0001 g

Serial No.: 1129361010

ID No.: UAE.WAS.002/2552

Max. Capacity: 220 g

Calibration Date: April 23, 2025

Condition As-Received: In Condition

Condition of Equipment:

Condition of This Result of Calibration:

1. Calibration Method: This instrument was calibrated by method UAE.CP.CAL.006 In-House Method based on UKAS Lab 14 : 2022

2. Reference Standards:

Reference Standard:	Model	Serial No.	Calibrated By	Certificate No.	Traceability	Due Date
Standard Weight Class E2 (OIML)	1 mg to 1 kg	B749109122	AMARC	25-009359	Mettler-Toledo	21-Jan-27
Standard Weight Class F1 (OIML)	1 mg to 200 g	11119512	AMARC	24-013840	Mettler-Toledo	04-Feb-26
Instrument	Model	Serial No.	Calibrated By	Certificate No.	Traceability	Due Date
Thermo-Hygro-Baro Meter	MHB-382SD	AK.46457	SUCCESS	SG-H-00997/67	Success Gateway	21-Nov-25
Thermo-Hygro-Baro Meter	MHB-382SD	AK.46457	TPA	25P795	TPA	25-Feb-26

3. This certification is traceable to SI Unit

4. This certification was certified only for the indtrument we calibrated

5. This result of calibration wae found accurate as show on date and place of calibration only.

6. Through the reference standard laboratory of AMARC 25-009359 Calibration 0152

Calibraton Result:




1. Repeatability of Reading:

Nominal Value (g)	Standard Deviation of Reading (g)
200*	0.000045

2. Eccentric or off-center loading

A mass of 100 g was placed and moved to various position on pan

The Balance reading obtained is given in the table.

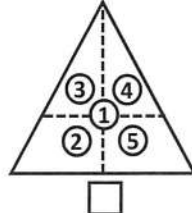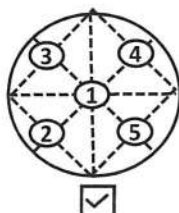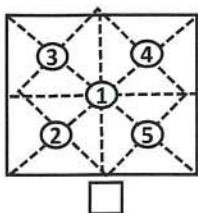

1 (g)	2 (g)	3 (g)	4 (g)	5 (g)	Maximum Difference (g)
100.0000	99.9996	99.9997	100.0003	100.0005	0.0005

เอกสารไม่ควบคุม

Certificate No.: 250422-1-BL002-25

Code No.: BL002-25

Page: 3 of 3

Equipment: Electronic Balance
Model: AB204-S/FACT
Serial No.: 1129361010
Max. Capacity: 220 g
Calibration Date: April 23, 2025

Manufacturer: Mettler Toledo
Readability: 0.0001 g
ID No.: UAE.WAS.002/2552

Calibration Result: (Continued)

Calibration Range: 0 - 200 g

Calibration Adjustment: Internal Calibration

3. Error of indication from nominal or conventional mass value:

Nominal Value (g)	Reference Value (g)	Indication (g)	Correction (g)	Uncertainty (\pm mg)	Coverage Factor <i>k</i>
Unload	0.0000000	0.0000	0.0000	0.10	2.05
0.01	0.0100025	0.0099	0.0001	0.10	2.05
0.05	0.0500056	0.0500	0.0000	0.10	2.05
0.1	0.1000012	0.0999	0.0001	0.10	2.05
0.5	0.5000133	0.5000	0.0000	0.10	2.05
1	1.0000105	1.0000	0.0000	0.10	2.05
10	10.000010	10.0000	0.0000	0.11	2.04
40	40.000076	40.0000	0.0000	0.14	2.00
50	50.000056	50.0000	0.0001	0.13	2.00
80	80.000107	80.0000	0.0001	0.18	2.00
100	100.000109	99.9999	0.0002	0.17	2.00
120	120.00015	119.9999	0.0003	0.21	2.00
150	150.000165	149.9998	0.0003	0.24	2.00
160	160.000175	159.9997	0.0005	0.26	2.00
200	200.000129	199.9998	0.0004	0.30	2.00

4. Effect of Tare test:

Tare Load (g)	Test Load (g)	Indication (g)	Correction (g)
100	20.000041	19.9999	0.0001
	40.000076	39.9998	0.0002
	60.000066	59.9997	0.0003
	80.000107	79.9999	0.0002
	100.000168	100.0004	-0.0003

Remark:

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor *k*, providing a level of confidence of approximately 95%.

เอกสารไม่ควบคุม

o---o-End-o---o

Calibration Certificate

Certificate No.: 2401718-001-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road,
Bangchack, Prakhonong, Bangkok 10260

Page 1 of 5

Equipment: pH Meter
Manufacturer: METTLER TOLEDO
Model: SevenEasy pH
Serial No.: 1231155210
ID No.: UAE.WAT.010/2553
Order No.: 2401718
Operation No.: 2401718-001
Date of Receipt: 27 February 2024
Date of Calibration: 11 March 2024

Calibrated by Mr.Manas Somsak
Specialist

Approved by 
(Mr.Pheraphat Tuanjit)
Manager, Division of Calibration Laboratory

Date of Issue: 12 March 2024

Responsible for the Technical Management Team

The uncertainties are for a confidence probability of approximately 95%.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65



Calibration Report

Certificate No.: 2401718-001-01

Equipment:

pH Meter

Resolution: 0.01 pH ; 1 mV

Manufacturer: METTLER TOLEDO

Model: SevenEasy pH

Serial No.: 1231155210

Type: Bench top

ID No.: UAE.WAT.010/2553

Date of Calibration: 11 March 2024

Page 2 of 5

Location: Chemical Calibration Laboratory, National Food Institute

Environment Condition: **Ambient Temperature:** (23.4 ± 1.5) °C **Relative Humidity:** (51 ± 3) %

Condition of Equipment: Good Condition

Condition of this Results of Calibration

1. Calibration Method W-CC-002 : In house method based on direct measurement by using standard voltage calibrator and certified reference material (CRM)

2. Reference Standards / Certified Reference Material

Instruments	Serial / ID No.	Manufacturer	Certificate No.	Due Date
2.1 DC Voltage Calibrator	2709007	Fluke	23E2003	14 June 2024
2.2 Digital Thermometer	2709007	Fluke	CC 660570-01	30 October 2024
2.3 Thermo-Hygro Meter	NFI.BTH 014/23	testo	CC 660353-01	3 April 2024
Certified Reference Material	Lot. No.	Manufacturer	Ref N	Expire Date
2.4 pH buffer 4.008 (Primary pH buffer Solution)	888842	CPAchem	PH216.L5	13 April 2025
2.5 pH buffer 6.865 (Primary pH buffer Solution)	888843	CPAchem	PH217.L5	13 April 2025
2.6 pH buffer 10.01 (Primary pH buffer Solution)	888844	CPAchem	PH220.L5	13 April 2024
2.7 pH buffer 7.00 (Standard pH buffer Solution)	C03109	HACH LANGE GmbH	S11M004	16 October 2025

3. This certification is traceable to The International System of Unit (SI Unit)

3.1 Instruments No. 2.1	through	NSC-TISI-TIS 17025 Laboratory Accreditation of Calibration No.0008
3.2 Instruments No. 2.2 and 2.3	through	NSC-TISI-TIS 17025 Laboratory Accreditation of Calibration No.0061
3.3 Certified Reference Material No. 2.4 to 2.6	traceable to	Primary measurement method- Harned cell using calibrated thermometer, barometer, and nanovoltmeter. The Standard Solution preparation and certified by CPAchem Ltd is accredited to ISO 17034 and ISO/IEC 17025
3.4 Certified Reference Material No. 2.7	traceable to	PTB Certificate Nr. PTB-PHOA-563/30504/23 and Certificate Nr. PTB-PHOB-555/30620/22 (PTB: Physikalisch-Technische Bundesanstalt, Braunschweig, Germany)

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

F-CS-012 Revision: 01 Date: 20-04-65

Calibration Report

Certificate No.: 2401718-001-01

Equipment: pH Meter **Resolution:** 0.01 pH ; 1 mV

Manufacturer: METTLER TOLEDO **Model:** SevenEasy pH

Serial No.: 1231155210 **Type:** Bench top

ID No.: UAE.WAT.010/2553

Date of Calibration: 11 March 2024

Page 3 of 5

Calibration Results:

1. Calibration of pH Meter (Manual Temperature Compensation at 25 °C)

(offset value before adjust: -0.4 mV)

Nominal pH	DC Voltage Standard (mV)	Average Indicator Reading		Uncertainty (±mV)	Coverage Factor (k)
		mV	pH		
0	414.121	414	0.00	0.58	2.00
2	295.814	296	2.00	0.58	2.00
4	177.464	178	4.00	0.58	2.00
6	59.160	59	6.00	0.58	2.00
7	0.001	0	7.00	0.58	2.00
8	-59.159	-59	8.00	0.58	2.00
10	-177.461	-177	10.00	0.58	2.00
12	-295.811	-296	12.00	0.58	2.00
14	-414.118	-414	14.00	0.58	2.00

2. Calibration of pH Meter with Electrode (Manual Temperature Compensation at 25 °C)

Equipment: pH Electrode **Type:** Combined Electrode

Manufacturer: METTLER TOLEDO **Model:** InLab Solids

Serial No.: 3065701 **ID.No.:** N/A

Performance of Electrode system (Three-Point Calibration at pH 4, 7 and 10)

Certified Value @25 °C (pH)	Average Indicator Reading		Relative Slope (%)	Uncertainty (± pH)	Coverage Factor (k)
	pH	mV			
4.008	4.01	188	-	0.0071	2.00
7.001	7.00	13	98.9	0.0086	2.00
10.010	10.01	-160	97.2	0.0085	2.00
6.865	6.87	21	-	0.0074	2.00

Calibration Report

Certificate No.: 2401718-001-01

Equipment: Digital Thermometer with RTD (pH Meter)

Resolution: 0.1 °C **Model:** SevenEasy pH
Serial No.: 1231155210 **ID No.:** UAE.WAT.010/2553
Manufacturer: METTLER TOLEDO

Date of Calibration: 11 March 2024

Page 4 of 5

Location: Chemical Calibration Laboratory, National Food Institute

Environment Condition:
Ambient Temperature 23 °C ± 1 °C
Relative Humidity 51 % ± 2 %

Condition of this results of Calibration:

1. Calibration Method :
 - In house method: W-TE-025 by comparison with standard thermometer.
 - The Calibration is determined by comparing with a known temperature from a standard resistance thermometer.
 - The temperature scale in use at this laboratory is the International Temperature scale of 1990 (ITS-90).

2. Reference Standard Instrument :

Instrument	Model	Serial No.	Certificate No.	Due Date	Through
HANDHELD THERMOMETER	1523	2118154	PSL-T 0877/66	06-Jun-24	TISTR
Platinum Resistance Thermometer (PRT)	5627A	877332			

Support Equipment : - Low Temperature Bath (ISOCAL-6), Model: Europa-6 Plus Basic, S/N: 341592/2

3. This certificate is traceable to International System of Units (SI Units).
4. This certificate was certified only for the instrument we calibrated.
5. This result of calibration was found accurate as shown on date and place of calibration only.

6. Condition of Calibrated item : Good

7. Result of Calibration :

☒

Without adjustment

☐

After adjustment



F-CS-012 Revision: 01 Date: 20-04-65



Calibration Report

Certificate No.: 2401718-001-01

Equipment: Digital Thermometer with RTD (pH Meter)

Resolution: 0.1 °C Model: SevenEasy pH

Serial No.: 1231155210 ID No.: UAE.WAT.010/2553

Manufacturer: METTLER TOLEDO

Date of Calibration: 11 March 2024

Page 5 of 5

Calibration point: 15.0, 25.0 and 35.0 °C

Calibration result:

- The probe was immersed in liquid bath or dry bath to a minimum depth of 100 mm.
- Description of probe, model : N/A S/N : N/A
- Dimension of probe : Diameter 4 mm., Length 120 mm.,
- Sheath material : Stainless Steel

UUC* Reading (°C)	Standard Temperature (°C)	Correction Value (°C)	Uncertainty ± (°C)
15.1	14.998	0.1	0.099
25.1	24.998	0.1	0.099
* 35.1	34.997	0.1	0.099

Note

- UUC* : Unit Under Calibration

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k= 2, providing a level of confidence of approximately 95 %.

----- End -----



F-CS-012 Revision: 01 Date: 20-04-65



Calibration Certificate

Certificate No.: 2402283-001-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO., LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road,
Bangchack, Prakhnong, Bangkok 10260

Page 1 of 4

Equipment: Electronic Balance
Manufacturer: METTLER TOLEDO
Model: XSR205DU
Serial No.: C009071872
ID No.: UAE.WAO.012/2563
Order No.: 2402283
Operation No.: 2402283-001
Date of Receipt: 2 April 2024
Date of Calibration: 2 April 2024

Calibrated by Mr.Jerawut Prapawuttipong
Scientist

Approved by



(Mr.Pheraphat Tuanjit)

Manager, Division of Calibration Laboratory

Date of Issue: 9 April 2024

Responsible for the Technical Management Team

The uncertainties are for a confidence probability of approximately 95%

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65



Calibration Report

Certificate No.: 2402283-001-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C009071872

ID No.: UAE.WAO.012/2563

Capacity: 220 g

Date of Calibration: 2 April 2024

Page 2 of 4

Environment Condition: Ambient Temperature: 24.5 ± 0.5 °C Relative Humidity: 47.5 ± 2.5 %

Place of Calibration: Laboratory, UNITED ANALYST AND ENGINEERING CONSULTANT CO., LTD.

Condition of Equipment: Good Condition

Condition of This Results of Calibration:

1. Calibration Method: NFI Method W-MA-001 In-House Method based on UKAS Lab 14 : 2019

2. Reference Standards:

Reference Standard	Model	Serial No.	Calibrated By	Certificate No.	Due Date
Standard Weight Class E2	1mg to 200g	B505567572	TCS	M2304053S	8 April 2024

Instrument	Model	Serial No.	Calibrated By	Certificate No.	Due Date
Thermo-Hygro Meter	608-H1	NFI.BTH 016/23	Quality Reborn	QR24-0343	9 February 2025

3. This certification is traceable to SI UNIT

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

Calibration Results:




1. Repeatability of Reading:

Nominal Value (g)	Standard Deviation of Reading (g)
40	0.0000052
80	0.0000063
100	0.000048
200	0.000053

2. Off-Center Error:

A mass of 100 g was placed and moved to various position on pan.

The balance reading obtained is given in the table.

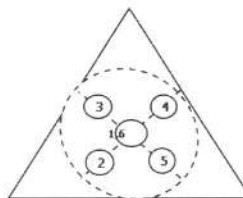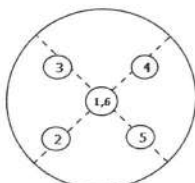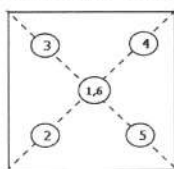

1 (g)	2 (g)	3 (g)	4 (g)	5 (g)	6 (g)	(Maximum Difference) (g)
100.0002	100.0001	100.0002	99.9999	100.0001	100.0001	0.0003



F-CS-012 Revision: 01 Date: 20-04-65

Calibration Report

Certificate No.: 2402283-001-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C009071872

ID No.: UAE.WAO.012/2563

Capacity: 220 g

Date of Calibration: 2 April 2024

Page 3 of 4

Calibration Results: (Continued)

Calibration Range: 0 - 80 g

Calibration Adjustment: Internal Calibration

3. Departure from Nominal Value: (Range: 0 - 80 g ; Resolution: 0.00001 g)

Nominal Value (g)	Standard Value (g)	Average Reading (g)	Correction (g)	Uncertainty (± g)	Coverage Factor k
Unload	0.000000	0.00000	0.00000	0.0000088	2.00
0.001	0.001003	0.00101	-0.00001	0.0000091	2.00
0.005	0.005003	0.00499	0.00001	0.0000094	2.00
0.01	0.010003	0.01000	0.00000	0.0000091	2.00
0.05	0.049996	0.05000	0.00000	0.0000098	2.00
0.1	0.100011	0.10000	0.00001	0.000011	2.00
0.5	0.500016	0.50001	0.00001	0.000014	2.00
1	1.000003	1.00002	-0.00002	0.000016	2.00
2	2.000023	2.00001	0.00001	0.000017	2.00
5	5.000017	5.00002	0.00000	0.000020	2.00
10	10.000009	10.00000	0.00001	0.000026	2.00
20	20.000031	20.00002	0.00001	0.000037	2.00
30	30.000040	30.00003	0.00001	0.000052	2.00
50	50.000028	50.00004	-0.00001	0.000068	2.00
80	80.000068	80.00005	0.00002	0.00011	2.00



Calibration Report

Certificate No.: 2402283-001-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C009071872

ID No.: UAE.WAO.012/2563

Capacity: 220 g

Date of Calibration: 2 April 2024

Page 4 of 4

Calibration Results: (Continued)

Calibration Range: 81 - 200 g

Calibration Adjustment: Internal Calibration

3. Departure from Nominal Value: (Range: 81 - 200 g ; Resolution: 0.0001 g)

Nominal Value (g)	Standard Value (g)	Average Reading (g)	Correction (g)	Uncertainty (± g)	Coverage Factor <i>k</i>
90	90.00010	90.0000	0.0001	0.00015	2.00
100	100.00006	100.0000	0.0001	0.00015	2.00
110	110.00007	110.0001	0.0000	0.00017	2.00
120	120.00009	120.0000	0.0001	0.00018	2.00
130	130.00010	130.0000	0.0001	0.00019	2.00
140	140.00014	140.0000	0.0001	0.00020	2.00
150	150.00009	150.0001	0.0000	0.00020	2.00
160	160.00010	160.0001	0.0000	0.00022	2.00
170	170.00012	170.0001	0.0000	0.00023	2.00
200	200.00016	200.0000	0.0002	0.00028	2.00

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k* , providing a level of confidence of approximately 95 %.

----- End -----

F-CS-012 Revision: 01 Date: 20-04-65



Calibration Certificate

Certificate No.: 2402283-002-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO., LTD.
Address: 3 SOI UDOMSUK 41, SUKHUMVIT ROAD,
Bangchack, Prakhonong, Bangkok 10260

Page 1 of 4

Equipment: Electronic Balance
Manufacturer: METTLER TOLEDO
Model: XSR205DU
Serial No.: C210685394
ID No.: UAE.WAO.010/2565
Order No.: 2402283
Operation No.: 2402283-002
Date of Receipt: 2 April 2024
Date of Calibration: 2 April 2024

Calibrated by Mr.Jerawut Prapawuttipong
Scientist

Approved by

(Mr.Pheraphat Tuanjit)

Manager, Division of Calibration Laboratory

Date of Issue: 9 April 2024

Responsible for the Technical Management Team

The uncertainties are for a confidence probability of approximately 95%

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65

Calibration Report

Certificate No.: 2402283-002-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C210685394

ID No.: UAE.WAO.010/2565

Capacity: 220 g

Date of Calibration: 2 April 2024

Page 2 of 4

Environment Condition: Ambient Temperature: 24.5 ± 0.5 °C Relative Humidity: 47.5 ± 2.5 %

Place of Calibration: Laboratory, UNITED ANALYST AND ENGINEERING CONSULTANT CO., LTD.

Condition of Equipment: Good Condition

Condition of This Results of Calibration:

1. Calibration Method: NFI Method W-MA-001 In-House Method based on UKAS Lab 14 : 2019

2. Reference Standards:

Reference Standard	Model	Serial No.	Calibrated By	Certificate No.	Due Date
Standard Weight Class E2	1mg to 200g	B505567572	TCS	M2304053S	8 April 2024

Instrument	Model	Serial No.	Calibrated By	Certificate No.	Due Date
Thermo-Hygro Meter	608-H1	NFI.BTH 016/23	Quality Reborn	QR24-0343	9 February 2025

3. This certification is traceable to SI UNIT

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

Calibration Results:




1. Repeatability of Reading:

Nominal Value (g)	Standard Deviation of Reading (g)
40	0.0000042
80	0.0000052
100	0.000048
200	0.000048

2. Off-Center Error:

A mass of 100 g was placed and moved to various position on pan.

The balance reading obtained is given in the table.

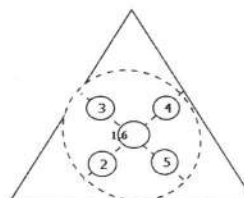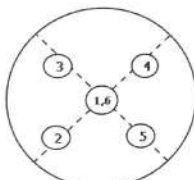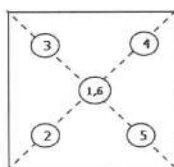

1 (g)	2 (g)	3 (g)	4 (g)	5 (g)	6 (g)	(Maximum Difference) (g)
100.0000	100.0001	99.9999	99.9999	100.0001	100.0000	0.0001




Calibration Report

Certificate No.: 2402283-002-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C210685394

ID No.: UAE.WAO.010/2565

Capacity: 220 g

Date of Calibration: 2 April 2024

Page 3 of 4

Calibration Results: (Continued)

Calibration Range: 0 - 80 g

Calibration Adjustment: Internal Calibration

3. Departure from Nominal Value: (Range: 0 - 80 g ; Resolution: 0.00001 g)

Nominal Value (g)	Standard Value (g)	Average Reading (g)	Correction (g)	Uncertainty (± g)	Coverage Factor k
Unload	0.000000	0.00000	0.00000	0.0000086	2.00
0.001	0.001003	0.00101	-0.00001	0.0000089	2.00
0.005	0.005003	0.00500	0.00000	0.0000092	2.00
0.01	0.010003	0.01000	0.00000	0.0000089	2.00
0.05	0.049996	0.05000	0.00000	0.0000096	2.00
0.1	0.100011	0.10000	0.00001	0.000011	2.00
0.5	0.500016	0.50001	0.00001	0.000014	2.00
1	1.000003	1.00002	-0.00002	0.000016	2.00
2	2.000023	2.00001	0.00001	0.000017	2.00
5	5.000017	5.00002	0.00000	0.000020	2.00
10	10.000009	10.00000	0.00001	0.000026	2.00
20	20.000031	20.00000	0.00003	0.000037	2.00
30	30.000040	30.00001	0.00003	0.000050	2.00
50	50.000028	50.00002	0.00001	0.000068	2.00
80	80.000068	80.00002	0.00005	0.00011	2.00



Calibration Report

Certificate No.: 2402283-002-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C210685394

ID No.: UAE.WAO.010/2565

Capacity: 220 g

Date of Calibration: 2 April 2024

Page 4 of 4

Calibration Results: (Continued)

Calibration Range: 81 - 200 g

Calibration Adjustment: Internal Calibration

3. Departure from Nominal Value: (Range: 81 - 200 g ; Resolution: 0.0001 g)

Nominal Value (g)	Standard Value (g)	Average Reading (g)	Correction (g)	Uncertainty (± g)	Coverage Factor k
90	90.00010	90.0001	0.0000	0.00015	2.00
100	100.00006	100.0001	0.0000	0.00015	2.00
110	110.00007	110.0001	0.0000	0.00016	2.00
120	120.00009	120.0000	0.0001	0.00017	2.00
130	130.00010	130.0000	0.0001	0.00019	2.00
140	140.00014	140.0000	0.0001	0.00020	2.00
150	150.00009	150.0001	0.0000	0.00020	2.00
160	160.00010	160.0001	0.0000	0.00022	2.00
170	170.00012	170.0001	0.0000	0.00023	2.00
200	200.00016	200.0002	0.0000	0.00028	2.00

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k , providing a level of confidence of approximately 95 %.

----- End -----

F-CS-012 Revision: 01 Date: 20-04-65



Calibration Certificate

Certificate No.: 2501844-001-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road,
Bangchack, Prakhonong, Bangkok 10260

Page 1 of 5

Equipment: pH Meter
Manufacturer: METTLER TOLEDO
Model: SevenEasy pH
Serial No.: 1231155210
ID No.: UAE.WAT.010/2553
Order No.: 2501844
Operation No.: 2501844-001
Date of Receipt: 24 February 2025
Date of Calibration: 24 February 2025

Calibrated by Mr.Manas Somsak
Specialist

Approved by 
(Mr.Pheraphat Tuanjit)
Manager, Division of Calibration Laboratory

Date of Issue: 27 February 2025

Responsible for the Technical Management Team

The uncertainties are for a confidence probability of approximately 95%.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65



Calibration Report

Certificate No.: 2501844-001-01

Equipment: pH Meter **Resolution:** 0.01 pH ; 1 mV

Manufacturer: METTLER TOLEDO **Model:** SevenEasy pH

Serial No.: 1231155210 **Type:** Bench top

ID No.: UAE.WAT.010/2553

Date of Calibration: 24 February 2025

Page 2 of 5

Location: Chemical Calibration Laboratory, National Food Institute

Environment Condition: **Ambient Temperature:** (23.4 ± 1.5) °C **Relative Humidity:** (54 ± 3) %

Condition of Equipment: Good Condition

Condition of this Results of Calibration

1. Calibration Method W-CC-002 : In house method based on direct measurement by using standard voltage calibrator and certified reference material (CRM)

2. Reference Standards / Certified Reference Material

Instruments	Serial / ID No.	Manufacturer	Certificate No.	Due Date
2.1 DC Voltage Calibrator	2709007	Fluke	24E1752	30 May 2025
2.2 Digital Thermometer	2709007	Fluke	2500376-002-01	29 October 2025
2.3 Thermo-Hygro Meter	NFI.BTH 013/23	testo	CC 670420-01	21 May 2025
Certified Reference Material	Lot. No.	Manufacturer	Ref N	Expire Date
2.4 pH buffer 4.008 (Primary pH buffer Solution)	1016435	CPAchem	PH216.L5	25 July 2026
2.5 pH buffer 6.865 (Primary pH buffer Solution)	949186	CPAchem	PH217.L5	30 November 2025
2.6 pH buffer 10.01 (Primary pH buffer Solution)	1016437	CPAchem	PH220.L5	25 July 2025
2.7 pH buffer 7.00 (Standard pH buffer Solution)	C03109	HACH LANGE GmbH	S11M004	16 October 2025

3. This certification is traceable to The International System of Unit (SI Unit)

3.1 Instruments No.2.1	through	NSC-TISI-TIS 17025 Laboratory Accreditation of Calibration No.0008
3.2 Instruments No.2.2 to 2.3	through	NSC-TISI-TIS 17025 Laboratory Accreditation of Calibration No.0061
3.3 Certified Reference Material No.2.4 to 2.6	traceable to	Primary measurement method- Harned cell using calibrated thermometer, barometer, and nanovoltmeter. The Standard Solution preparation and certified by CPAchem Ltd is accredited to ISO 17034 and ISO/IEC 17025
3.4 Certified Reference Material No.2.7	traceable to	PTB Certificate Nr. PTB-PHOA-563/30504/23 and Certificate Nr. PTB-PHOB-555/30620/22 (PTB: Physikalisch-Technische Bundesanstalt, Braunschweig, Germany)

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.



F-CS-012 Revision: 01 Date: 20-04-65



Calibration Report

Certificate No.: 2501844-001-01

Equipment: pH Meter **Resolution:** 0.01 pH ; 1 mV

Manufacturer: METTLER TOLEDO **Model:** SevenEasy pH

Serial No.: 1231155210 **Type:** Bench top

ID No.: UAE.WAT.010/2553

Date of Calibration: 24 February 2025

Page 3 of 5

Calibration Results:

1. Calibration of pH Meter (Manual Temperature Compensation at 25 °C)

Nominal pH	DC Voltage Standard (mV)	Average Indicator Reading		Uncertainty (±mV)	Coverage Factor (k)
		mV	pH		
0	414.122	414	-0.01	0.58	2.00
2	295.815	296	1.99	0.58	2.00
4	177.463	178	4.00	0.58	2.00
6	59.160	59	6.00	0.58	2.00
7	0.001	0	7.00	0.58	2.00
8	-59.159	-59	8.00	0.58	2.00
10	-177.462	-177	10.00	0.58	2.00
12	-295.813	-296	12.00	0.58	2.00
14	-414.121	-414	14.00	0.58	2.00

2. Calibration of pH Meter with Electrode (Manual Temperature Compensation at 25 °C)

Equipment: pH Electrode **Type:** Combined Electrode

Manufacturer: METTLER TOLEDO **Model:** InLab Solids

Serial No.: 3065701 **ID.No.:** N/A

Performance of Electrode system (Three-Point Calibration at pH 4, 7 and 10)

Certified Value @25 °C (pH)	Average Indicator Reading		Relative Slope (%)	Uncertainty (± pH)	Coverage Factor (k)
	pH	mV			
4.008	4.00	165	-	0.0071	2.00
7.001	7.00	-8	97.5	0.0086	2.00
10.010	10.01	-178	95.5	0.0083	2.00
6.876	6.88	0	-	0.0071	2.00



Calibration Report

Certificate No.: 2501844-001-01

Equipment: Digital Thermometer with RTD (pH Meter)

Resolution: 0.1 °C **Model:** SevenEasy pH
Serial No.: 1231155210 **ID No.:** UAE.WAT.010/2553
Manufacturer: METTLER TOLEDO

Date of Calibration: 24 February 2025

Page 4 of 5

Location: Chemical Calibration Laboratory, National Food Institute

Environment Condition:
Ambient Temperature 23.4 °C ± 1.0 °C
Relative Humidity 55.1 % ± 1.7 %

Condition of this results of Calibration:

1. Calibration Method :
 - In house method: W-TE-025 by comparison with standard thermometer.
 - The Calibration is determined by comparing with a known temperature from a standard resistance thermometer.
 - The temperature scale in use at this laboratory is the International Temperature scale of 1990 (ITS-90).

2. Reference Standard Instrument :

Instrument	Model	Serial No.	Certificate No.	Due Date	Through
HANDHELD THERMOMETER	1523	2118154	PSL-T 0815/67	24-Jun-25	TISTR
Platinum Resistance Thermometer (PRT)	5627A	877332			

Support Equipment : - Low Temperature Bath (AMETEK RTC-187), Model: RTC-187C , S/N: 670930-00018

3. This certificate is traceable to International System of Units (SI Units).
4. This certificate was certified only for the instrument we calibrated.
5. This result of calibration was found accurate as shown on date and place of calibration only.

6. Condition of Calibrated item : Good

7. Result of Calibration : ☒ Without adjustment ☐ After adjustment




Calibration Report

Certificate No.: 2501844-001-01

Equipment: Digital Thermometer with RTD (pH Meter)

Resolution: 0.1 °C Model: SevenEasy pH

Serial No.: 1231155210 ID No.: UAE.WAT.010/2553

Manufacturer: METTLER TOLEDO

Date of Calibration: 24 February 2025

Page 5 of 5

Calibration point: 20.0, 25.0 and 30.0 °C

Calibration result:

- The probe was immersed in liquid bath or dry bath to a minimum depth of 120 mm.
- Description of probe, model : N/A S/N : N/A
- Dimension of probe : Diameter 4 mm., Length 120 mm.,
- Sheath material : Stainless Steel

UUC* Reading (°C)	Standard Temperature (°C)	Correction Value (°C)	Uncertainty ± (°C)
20.1	20.001	0.1	0.099
25.1	25.002	0.1	0.099
30.1	30.003	0.1	0.099

Note

- UUC* : Unit Under Calibration

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k= 2, providing a level of confidence of approximately 95 %.

----- End -----



F-CS-012 Revision: 01 Date: 20-04-65


Calibration Certificate

Certificate No.: 2502226-001-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road,
Bangchack, Prakhnong, Bangkok 10260

Page 1 of 4

Equipment: Electronic Balance
Manufacturer: METTLER TOLEDO
Model: XSR205DU
Serial No.: C009071872
ID No.: UAE.WAO.012/2563
Order No.: 2502226
Operation No.: 2502226-001
Date of Receipt: 19 March 2025
Date of Calibration: 20 March 2025

Calibrated by Mr.Yothin Charoensuk
Scientist

Approved by 
(Mr.Pheraphat Tuanjit)
Manager, Division of Calibration Laboratory
Responsible for the Technical Management Team

Date of Issue: 25 March 2025

The uncertainties are for a confidence probability of approximately 95%

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65

Calibration Report

Certificate No.: 2502226-001-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C009071872

ID No.: UAE.WAO.012/2563

Capacity: 82 g / 220 g

Date of Calibration: 20 March 2025

Page 2 of 4

Environment Condition: Ambient Temperature: 21.2 ± 0.6 °C Relative Humidity: 48 ± 3.5 %

Place of Calibration: 208 Balance Room, UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.

Condition of Equipment: Good Condition

Condition of This Results of Calibration:

1. Calibration Method: NFI Method W-MA-001 In-House Method based on UKAS Lab 14 : 2019

2. Reference Standards:

Reference Standard	Model	Serial No.	Calibrated By	Certificate No.	Due Date
Standard Weight Class E2	1mg to 200g	B505567572	TCS	M2404100S	19 April 2025

Instrument	Model	Serial No.	Calibrated By	Certificate No.	Due Date
Thermo-Hygro Meter	608-H1	NFI.BTH 017/23	Quality Reborn	QR25-0542	10 February 2026

3. This certification is traceable to SI UNIT

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

Calibration Results:




1. Repeatability of Reading:

Nominal Value (g)	Standard Deviation of Reading (g)
40	0.0000052
80	0.0000042
100	0.0000000
200	0.0000000

2. Off-Center Error:

A mass of 100 g was placed and moved to various position on pan.

The balance reading obtained is given in the table.

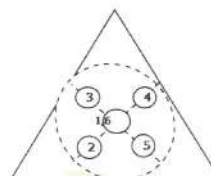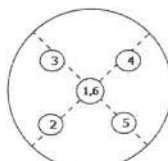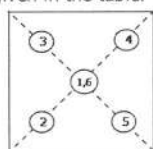

1	2	3	4	5	6	(Maximum Difference)
(g)	(g)	(g)	(g)	(g)	(g)	(g)
100.0001	100.0001	100.0001	100.0001	100.0001	100.0002	0.0001

for N. mityambut

F-CS-012 Revision: 01 Date: 20-04-65

Calibration Report

Certificate No.: 2502226-001-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C009071872

ID No.: UAE.WAO.012/2563

Capacity: 82 g / 220 g

Date of Calibration: 20 March 2025

Page 3 of 4

Calibration Results: (Continued)

Calibration Range: 0-80 g

Calibration Adjustment: Internal Calibration

3. Departure from Nominal Value: (Range: 0 - 82 g ; Resolution: 0.00001 g)

Nominal Value (g)	Standard Value (g)	Average Reading (g)	Correction (g)	Uncertainty (± g)	Coverage Factor k
Unload	0.000000	0.00000	0.00000	0.0000089	2.00
0.001	0.001003	0.00100	0.00000	0.0000092	2.00
0.005	0.005002	0.00500	0.00000	0.0000094	2.00
0.01	0.010003	0.01000	0.00000	0.0000091	2.00
0.05	0.049996	0.05000	0.00000	0.0000098	2.00
0.1	0.100011	0.10000	0.00001	0.000011	2.00
0.5	0.500016	0.50000	0.00002	0.000014	2.00
1	1.000003	1.00001	-0.00001	0.000016	2.00
2	2.000023	2.00005	-0.00003	0.000017	2.00
5	5.000015	5.00005	-0.00003	0.000021	2.00
10	10.000009	10.00005	-0.00004	0.000026	2.00
20	20.000030	20.00012	-0.00009	0.000037	2.00
30	30.000039	30.00012	-0.00008	0.000050	2.00
50	50.000028	50.00014	-0.00011	0.000068	2.00
80	80.000067	80.00020	-0.00013	0.00011	2.00

Calibration Report

Certificate No.: 2502226-001-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C009071872

ID No.: UAE.WAO.012/2563

Capacity: 82 g / 220 g

Date of Calibration: 20 March 2025

Page 4 of 4

Calibration Results: (Continued)

Calibration Range: >80-200 g

Calibration Adjustment: Internal Calibration

3. Departure from Nominal Value: (Range: >80 - 200 g ; Resolution: 0.0001 g)

Nominal Value (g)	Standard Value (g)	Average Reading (g)	Correction (g)	Uncertainty (± g)	Coverage Factor <i>k</i>
90	90.00010	90.0002	-0.0001	0.00015	2.00
100	100.00006	100.0001	0.0000	0.00016	2.00
110	110.00007	110.0001	0.0000	0.00017	2.00
120	120.00009	120.0002	-0.0001	0.00018	2.00
130	130.00010	130.0002	-0.0001	0.00019	2.00
140	140.00013	140.0002	-0.0001	0.00019	2.00
150	150.00009	150.0002	-0.0001	0.00021	2.00
160	160.00010	160.0002	-0.0001	0.00022	2.00
170	170.00012	170.0002	-0.0001	0.00023	2.00
200	200.00013	200.0002	-0.0001	0.00028	2.00

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k* , providing a level of confidence of approximately 95 %.

----- End -----

F-CS-012 Revision: 01 Date: 20-04-65

for N. ingudat


Calibration Certificate

Certificate No.: 2502226-002-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road,
Bangchack, Prakhnong, Bangkok 10260

Page 1 of 4

Equipment: Electronic Balance
Manufacturer: METTLER TOLEDO
Model: XSR205DU
Serial No.: C210685394
ID No.: UAE.WAO.010/2565
Order No.: 2502226
Operation No.: 2502226-002
Date of Receipt: 19 March 2025
Date of Calibration: 20 March 2025

Calibrated by Mr.Yothin Charoensuk
Scientist

Approved by 
(Mr.Pheraphat Tuanjit)
Manager, Division of Calibration Laboratory
Responsible for the Technical Management Team

Date of Issue: 25 March 2025

The uncertainties are for a confidence probability of approximately 95%

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65

Calibration Report

Certificate No.: 2502226-002-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C210685394

ID No.: UAE.WAO.010/2565

Capacity: 82 g / 220 g

Date of Calibration: 20 March 2025

Page 2 of 4

Environment Condition: Ambient Temperature: 21.2 ± 0.6 °C Relative Humidity: 48 ± 3.5 %

Place of Calibration: 208 Balance Room, UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.

Condition of Equipment: Good Condition

Condition of This Results of Calibration:

1. Calibration Method: NFI Method W-MA-001 In-House Method based on UKAS Lab 14 : 2019

2. Reference Standards:

Reference Standard	Model	Serial No.	Calibrated By	Certificate No.	Due Date
Standard Weight Class E2	1mg to 200g	B505567572	TCS	M2404100S	19 April 2025

Instrument	Model	Serial No.	Calibrated By	Certificate No.	Due Date
Thermo-Hygro Meter	608-H1	NFI.BTH 017/23	Quality Reborn	QR25-0542	10 February 2026

3. This certification is traceable to SI UNIT

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

Calibration Results:




1. Repeatability of Reading:

Nominal Value (g)	Standard Deviation of Reading (g)
40	0.0000042
80	0.0000042
100	0.0000000
200	0.0000000

2. Off-Center Error:

A mass of 100 g was placed and moved to various position on pan.

The balance reading obtained is given in the table.

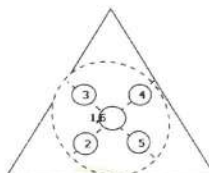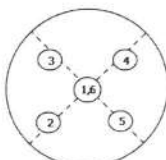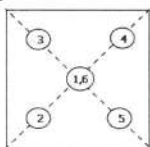

1	2	3	4	5	6	(Maximum Difference)
(g)	(g)	(g)	(g)	(g)	(g)	(g)
100.0001	100.0001	100.0001	100.0001	100.0001	100.0001	0.0000

Calibration Report

Certificate No.: 2502226-002-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C210685394

ID No.: UAE.WAO.010/2565

Capacity: 82 g / 220 g

Date of Calibration: 20 March 2025

Page 3 of 4

Calibration Results: (Continued)

Calibration Range: 0-80 g

Calibration Adjustment: Internal Calibration

3. Departure from Nominal Value: (Range: 0 - 82 g ; Resolution: 0.00001 g)

Nominal Value (g)	Standard Value (g)	Average Reading (g)	Correction (g)	Uncertainty (± g)	Coverage Factor k
Unload	0.000000	0.00000	0.00000	0.0000087	2.00
0.001	0.001003	0.00100	0.00000	0.0000090	2.00
0.005	0.005002	0.00501	-0.00001	0.0000092	2.00
0.01	0.010003	0.01002	-0.00002	0.0000089	2.00
0.05	0.049996	0.05001	-0.00001	0.0000096	2.00
0.1	0.100011	0.10002	-0.00001	0.000011	2.00
0.5	0.500016	0.50004	-0.00002	0.000014	2.00
1	1.000003	1.00005	-0.00005	0.000016	2.00
2	2.000023	2.00006	-0.00004	0.000017	2.00
5	5.000015	5.00006	-0.00005	0.000020	2.00
10	10.000009	10.00005	-0.00004	0.000026	2.00
20	20.000030	20.00007	-0.00004	0.000037	2.00
30	30.000039	30.00009	-0.00005	0.000050	2.00
50	50.000028	50.00008	-0.00005	0.000068	2.00
80	80.000067	80.00013	-0.00006	0.00011	2.00



Calibration Report

Certificate No.: 2502226-002-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C210685394

ID No.: UAE.WAO.010/2565

Capacity: 82 g / 220 g

Date of Calibration: 20 March 2025

Page 4 of 4

Calibration Results: (Continued)

Calibration Range: >80-200 g

Calibration Adjustment: Internal Calibration

3. Departure from Nominal Value: (Range: >80 - 200 g ; Resolution: 0.0001 g)

Nominal Value (g)	Standard Value (g)	Average Reading (g)	Correction (g)	Uncertainty (± g)	Coverage Factor <i>k</i>
90	90.00010	90.0002	-0.0001	0.00015	2.00
100	100.00006	100.0001	0.0000	0.00016	2.00
110	110.00007	110.0002	-0.0001	0.00017	2.00
120	120.00009	120.0002	-0.0001	0.00018	2.00
130	130.00010	130.0002	-0.0001	0.00019	2.00
140	140.00013	140.0002	-0.0001	0.00019	2.00
150	150.00009	150.0002	-0.0001	0.00021	2.00
160	160.00010	160.0002	-0.0001	0.00022	2.00
170	170.00012	170.0002	-0.0001	0.00023	2.00
200	200.00013	200.0002	-0.0001	0.00028	2.00

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k* , providing a level of confidence of approximately 95 %.

----- End -----

for N. important

F-CS-012 Revision: 01 Date: 20-04-65



Calibration Certificate

Certificate No.: 2502229-002-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road,
Bangchack, Prakhanong, Bangkok 10260

Page 1 of 3

Equipment: CHAMBER (Incubator)

Manufacturer: MEMMERT

Model: IPP260

Serial No.: V616.0066

ID No.: UAE.MIC.032/2559


Order No.: 2502229

Operation No.: 2502229-002

Date of Receipt: 19 March 2025

Date of Calibration: 19 March 2025

Calibrated by Mr.Yothin Charoensuk
Scientist

Approved by 
(Mr.Pheraphat Tuanjit) (for)
Manager, Division of Calibration Laboratory
Responsible for the Technical Management Team

Date of Issue: 25 March 2025

The uncertainties are for a confidence probability of approximately 95 %.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.



Calibration Report

Certificate No.: 2502229-002-01

Equipment: CHAMBER (Incubator)

Model: IPP260 Serial No.: V616.0066

Resolution: 0.1 °C ID No.: UAE.MIC.032/2559

Manufacturer: MEMMERT

Date of Calibration: 19 March 2025

Page 2 of 3

Location: 302, UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.

Environment Condition:
Ambient Temperature (21.7 ± 1) °C
Relative Humidity (59 ± 1) %
Line Voltage (223 ± 3) Volt

Condition of this results of Calibration:

- This instrument was calibrated by insert 9 standard thermometer into its chamber and calibration according to W-TE-014 Based on TLAS G-20-1/02-08 (E): Guidelines for Calibration and Checks of Temperature Controlled Enclosures.
- The temperature scale used was based on ITS - 90.
- All data show below were final values and the initial data may be obtained upon request.

2. Reference Standard Instrument :

Instrument	Model	Serial No./ID No.	Certificate No.	Due Date	Through
Digital Thermometer with sensor	34972A	MY57003188	TE 670486-01	8 June 2025	NATIONAL FOOD INSTITUTE
	RTD	CH#201-209/ RTD#201-209			

- This certificate is traceable to International System of Units (SI Units).
- This certificate was certified only for the instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.
- Condition of Calibrated item : Good

UUC Description :

Time of Record 1 Hour 9 Minute At 25.0 and 36.0 °C
Fresh air Damper ☐ Open Position ☐
☒ Close Fan ☐
☐ Not Available

7. Result of Calibration : ☒ Without adjustment ☐ After adjustment

P. Jongsakul
25 March 2025



Calibration Report

Certificate No.: 2502229-002-01

Equipment: CHAMBER (Incubator)

Model: IPP260 Serial No.: V616.0066

Resolution: 0.1 °C ID No.: UAE.MIC.032/2559

Manufacturer: MEMMERT

Date of Calibration: 19 March 2025

Page 3 of 3

Calibration point: 25.0 and 36.0 °C

Calibration result:

Calibration Condition	Temperature (°C)	Relative Humidity (%)	Line Voltage (Volt)
MIN	21.3	58	220.0
MAX	22.0	60	225.0

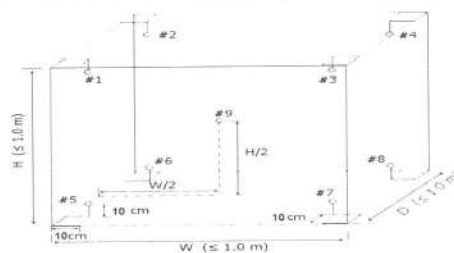


Table1 : Reporting of Temperature

Calibration point (°C)	Measured Temperature (°C) @ Sensor No. (Sensor No.9 is REF)									Uncertainty ± (°C)
	# 1	# 2	# 3	# 4	# 5	# 6	# 7	# 8	# 9	
25.0	25.19	25.16	25.22	25.17	24.85	24.91	24.78	24.85	24.97	0.29
36.0	34.57	34.74	35.13	35.29	36.32	36.16	36.20	36.34	35.73	0.63

Table 2 : Reporting of Characterization Result

UUC* Setting (°C)	UUC* Reading (°C)			Stability ± (°C)	Uniformity (°C)	Overall Variation (°C)
	MIN	MAX	Average			
25.0	24.9	25.0	25.0	0.088	0.25	0.61
36.0	35.9	36.0	36.0	0.44	1.2	2.3

Note The quoted uncertainty include " Stability " and " Loading effect (20% of Temp Uniformity) "

UUC* = Unit Under Calibration

Stability = One-half of the greatest maximum difference of measured temperatures at any one sensors, for at least half an hour after reaching steady state.

Uniformity = The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time.

Overall Variation = The difference of the maximum and minimum measured temperatures throughout observation time.

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor $k=2$, providing a level of confidence of approximately 95 %.

----- End -----

P. Jenghachit
25 March 2025



Calibration Certificate

Certificate No.: 2502229-006-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road,
Bangchack, Prakhonong, Bangkok 10260

Page 1 of 3

Equipment: CHAMBER (Incubator)
Manufacturer: BINDER
Model: KB 400
Serial No.: 20200000015535
ID No.: UAE.MIC.018/2564
Order No.: 2502229
Operation No.: 2502229-006
Date of Receipt: 19 March 2025
Date of Calibration: 19 March 2025

Calibrated by Mr.Jerawut Prapawuttipong
Scientist

Approved by


(Mr.Pheraphat Tuanjit) (for)

Manager, Division of Calibration Laboratory

Date of Issue: 25 March 2025

Responsible for the Technical Management Team

The uncertainties are for a confidence probability of approximately 95 %.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65



Calibration Report

Certificate No.: 2502229-006-01
Equipment: CHAMBER (Incubator)
Model: KB 400 Serial No.: 20200000015535
Resolution: 0.1 °C ID No.: UAE.MIC.018/2564
Manufacturer: BINDER
Date of Calibration: 19 March 2025

Page 2 of 3

Location: LABORATORY, UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Environment Condition: Ambient Temperature (18 ± 1) °C
Relative Humidity (50 ± 5) %
Line Voltage (223 ± 3) Volt

Condition of this results of Calibration:

- This instrument was calibrated by insert 13 standard thermometer into its chamber and calibration according to W-TE-014 Based on TLAS G-20-1/02-08 (E): Guidelines for Calibration and Checks of Temperature Controlled Enclosures.
- The temperature scale used was based on ITS - 90.
- All data show below were final values and the initial data may be obtained upon request.

2. Reference Standard Instrument :

Instrument	Model	Serial No./ID No.	Certificate No.	Due Date	Through
Digital Thermometer with sensor	34972A	MY49016851	TE 670477-01	4 May 2025	NATIONAL FOOD INSTITUTE
	RTD	CH#201-303 / RTD#201-303			

- This certificate is traceable to International System of Units (SI Units).
- This certificate was certified only for the instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.
- Condition of Calibrated item : Good

UUC Description :

Time of Record 1 Hour 9 Minute At 35.0 °C
Fresh air Damper ☐ Open Position ☐
☒ Close Fan ☐
☐ Not Available

7. Result of Calibration : ☒ Without adjustment ☐ After adjustment

P. Jongsakul
25 March 2025



Calibration Report

Certificate No.: 2502229-006-01
Equipment: CHAMBER (Incubator)
Model: KB 400 **Serial No.:** 20200000015535
Resolution: 0.1 °C **ID No.:** UAE.MIC.018/2564
Manufacturer: BINDER

Date of Calibration: 19 March 2025

Page 3 of 3

Calibration point: 35.0 °C

Calibration result:

Calibration Condition	Temperature (°C)	Relative Humidity (%)	Line Voltage (Volt)
MIN	17.1	45	220.0
MAX	18.1	55	225.0

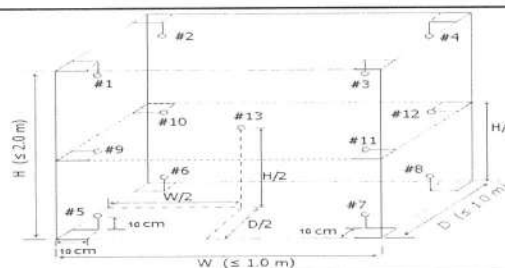


Table1 : Reporting of Temperature

Calibration point (°C)	Measured Temperature (°C) @ Sensor No. (Sensor No.13 is REF)													Uncertainty ± (°C)
	# 1	# 2	# 3	# 4	# 5	# 6	# 7	# 8	# 9	# 10	# 11	# 12	# 13	
35.0	34.98	35.17	34.99	34.92	35.18	35.01	35.00	35.13	35.00	34.96	35.02	35.17	35.04	0.27

Table 2 : Reporting of Characterization Result

UUC* Setting (°C)	UUC* Reading (°C)			Temperature Stability ± (°C)	Temperature Uniformity (°C)	Overall Variation (°C)
	MIN	MAX	Average			
35.0	35.0	35.0	35.0	0.029	0.15	0.30

Note The quoted uncertainty include " Stability " and " Loading effect (20% of Temp Uniformity) "

UUC* = Unit Under Calibration

Stability = One-half of the greatest maximum difference of measured temperatures at any one sensors, for at least half an hour after reaching steady state.

Uniformity = The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time.

Overall Variation = The difference of the maximum and minimum measured temperatures throughout observation time.

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k= 2, providing a level of confidence of approximately 95 %.

----- End -----

F-CS-012 Revision: 01 Date: 20-04-65

P. Jongsakul

25 March 2025

เอกสารไม่ควบคุม



Agilent CrossLab Start Up Services

Agilent 7890 Gas Chromatograph

Preventive Maintenance Checklist

Agilent Preventive Maintenance provides factory recommended service for your analytical instruments to assure reliable operation and the accuracy of your results.

Delivered by highly trained and certified service engineers using genuine Agilent parts and supplies, Agilent Preventive Maintenance provides everything you need to reduce unplanned downtime and keep your systems operating at their peak. This checklist will be completed at the end of the service and provided to you as a record of the preventive maintenance activities.

Introduction

Customer Information

- Customers should provide all necessary operating supplies upon request of the engineer.
- A customer representative should be available to the engineer while performing the preventive maintenance procedures.
- Any parts, not included in the Parts Lists section of this document, are not part of the recommended Preventive Maintenance service, nor are they included in the price of this service.
- If a system requires the use of extra or special procedures and/or parts for the maintenance service, then these must be ordered separately and charged as a repair, which may incur additional costs.

Important Customer Web Links

- For more information about **Agilent Technologies services**, please visit our website using the following URL: <http://www.agilent.com/en-us/products/crosslab-instrument-services/service-repair>
- The **Agilent Community** is an excellent place to get answers, collaborate with others about applications and Agilent products, and find in-depth documents and videos relevant to Agilent technologies. Visit <https://community.agilent.com/welcome>.
- To access **Agilent University**, visit <http://www.agilent.com/crosslab/university/> to learn about training options, which include online, classroom and onsite delivery. A training specialist can work directly with you to help determine your best options.
- A useful **Agilent Resource Center** web page is available, which includes short videos on maintenance, quick lists of consumables for new instruments, and other valuable information. Check out the Resource Page here: <https://www.agilent.com/en-us/agilentresources>.
- Need technical support, FAQs, supplies? – visit our **Support Home page** <http://www.agilent.com/search/support>.
- **Videos** about specific preparation requirements for your instrument can be found by searching the **Agilent YouTube** channel at <https://www.youtube.com/user/agilent>.
- **7890B Manuals** are also available on Agilent.com:
 - **Safety**
https://www.agilent.com/cs/library/usermanuals/public/7890B_Safety.pdf
 - **Installation and First Startup**
https://www.agilent.com/cs/library/usermanuals/Public/7890B_Installation.pdf
 - **Operation Manual**
https://www.agilent.com/cs/library/usermanuals/Public/7890B_Operation.pdf
 - **Maintaining Your GC**
https://www.agilent.com/cs/library/usermanuals/public/G3430-90052%207890B_Maintaining%20Guide.pdf

Service Engineer's Responsibilities

- Contact the customer and ensure that all necessary supplies are available before the preventive maintenance visit.
- Only select those pages that relate to the system or module being serviced.
- Complete empty fields with the relevant information.
- Complete the relevant checkboxes in the checklist using either a "X" or tick mark "✓".
- Check **"Section not applicable"** check boxes to indicate services/tasks not delivered, as appropriate.
- Complete the Preventive Maintenance service in the order of the tasks listed.
- Complete the Service Review section together with the customer.
- Complete the fields for page numbers at the foot of each selected page
- Complete the total number of pages field in the Service Completion section
- ***Ask the customer to sign the Service Completion section including the customer's and your signature.***

Additional Instruction Notes

- Check for any active service notes for this unit. If there are any applicable "Safety" or "Modification Recommended" Service notes, plan to implement the changes on this unit before doing any qualification service.
- Do not implement firmware updates, unless you get approval from the customer and are sure that they are compatible with the instrument control software.

System Information

- ☒ Check this box if an instrument configuration report is attached instead of completing the table below.

Instrument System Name and ID	UAE.TOX.007
Instrument System Site and Location	Laboratory

List System Component Product Numbers	List the Serial Numbers of each Component
1. G3440A	CN11021007
2. G4513A	CN20030059
3. G4514A	CN20020060
4.	
5.	
6.	
7.	
8.	
9.	
10.	

Preparation

- ☒ Discuss any specific issues with the customer before starting.
- ☒ Review the instrument logbook for recorded problems and comments.
- ☒ Save instrument control settings before starting the procedure.
- ☒ Perform a general inspection of the system for cleanliness.
- ☒ Check for proper installation of parts, assemblies, sensors etc.
- ☒ Check system for required installation of components, settings as defined by current Service Notes.
- ☒ Check for required firmware updates and verify with customers if they would like them installed.
- ☒ Before starting the following procedures, record the Detector Signal Output(s) in the results table. If the GC is turned OFF or in a service mode, comparing the detector outputs before and after the service is not possible.

Preventive Maintenance Procedure

Clean and inspect GC

- ☒ Unplug power cord from the power source.
- ☒ Open GC covers and vacuum/remove any dust/debris. Pay particular attention to cooling fans.
- ☒ Inspect internal connectors for proper contact and placement.
- ☒ Reconnect Power to the GC. Power the GC on and verify the power on self-test passed.
- ☒ Verify oven motor spins freely and turns on with the oven door closed; off when the door is opened.
- ☒ Verify operation of all other fans - the inlet and EPC cooling fans.
- ☒ Verify oven intake/outlet flap assembly is operating smoothly while heating and cooling the oven

Inlet and detector consumable replacement

- ☒ For the inlets installed, perform inlet maintenance as defined in the 7890 manual – “Maintaining Your GC” - for the inlet(s) installed.
- ☒ Replace the split vent trap cartridge filter on units with these inlets: Split/Splitless Capillary (SSL), Multi-Mode Inlet (MMI), Programmed Temperature Vaporizer (PTV), Volatiles Interface (VI).
- ☒ If the inlet system is used in Split Mode with viscous samples, inspect and clean the split vent tube on the inlet and flush or replace the tubing between the inlet and the split vent trap.
- ☒ If the GC includes a Flame Ionization Detector (FID), replace the jet. If the ignitor shows any buildup of sample or corrosion, replace the ignitor. Examine the FID collector and castle assemblies for contamination – clean as necessary.

Zero Sensors and Leak test

- ☒ Zero all pressure sensors per the procedure in the 7890 “Advanced User Guide”.
- ☒ Perform inlet pressure decay test(s) as defined in the 7890 “Troubleshooting Manual”.
If the PM is done in preparation for an Operational Qualification, then the pressure decay test defined within that protocol can be used for the PM.
- ☒ Record if test passed or failed in the results table.

ALS Maintenance

- ☐ **Section NOT applicable**
- ☒ Check all cabling and configuration settings between GC, tray, and injectors.
- ☒ Vacuum or remove any dust, especially around fans.
- ☒ Check operation of all fans.
- ☒ Check syringe for smooth plunger operation.
- ☒ Check for smooth operation of the needle support – clean if necessary

Restore Instrument

- ☒ Restore the normal operating conditions or customer method using the Data System.
- ☒ Purge the system with carrier flow for 15 minutes
- ☒ Bake out the system, then restore the normal operating conditions
- ☒ After equilibration, check and record the post PM detector signal output values.
Results should be similar or lower than the detector outputs recorded prior to PM.
- ☒ Perform a chemical checkout. If this is a routine PM, inject the customer's sample using the ALS if applicable. This will act as a final checkout of both the ALS and the GC.

Note: If the PM Service is performed prior to a qualification service, then use the qualification procedure as a guide for final instrument set up and checkout.

Signature Page

Service Review

- ☒ Attach available reports/printouts of all tests to this documentation.
- ☒ Record the Preventive Maintenance service activity in the customer's records/logbook.
- ☒ Update/reset instrument maintenance counters as appropriate.
- ☒ Affix the PM sticker to the system or instrument logbook based on the customer's request.
- ☒ Complete the Service Engineer Comments section if there are additional comments.
- ☒ Review with the customer this service, parts replaced, and test results obtained.
- ☐ If the instrument firmware was updated, record the details of the change in the Service Engineer's Comments box or if necessary, in the customer's IQ records.
- ☐ Supply the customer with a copy of the Smart Alerts flyer.
- ☐ Describe Smart Alerts to the customer.
- ☐ Install Smart Alerts if requested.

7890 GC Test Results Table

Detector Signal Outputs	Before PM Service	After PM Service
Front detector output	313.5	313.5
Back detector output	24.5	19.3
AUX detector output	n/a	n/a
Pressure decay test	Expected test result	Actual test result
Front inlet pressure decay test	Pass	Pass
Back inlet pressure decay test	Pass	n/a

7890 Parts List Table

The following kits are recommended for capillary and purged packed inlets. If this is a general PM and the customer has a preferred set of consumables, you may use the customer's consumables.

Part description	Part number	Product or model# where used	Quantity consumed
SSL Capillary Inlet PM kit, Splitless	5188-6497	7890A/B	1
SSL Capillary Inlet PM kit, split	5188-6496	7890A/B	-
SSL Capillary Ultra Inert Inlet Gold Seal with Washer	5190-6144	7890A/B	-
SSL Capillary Ultra Inert Inlet Splitless Liner - Single taper with Glass Wool	5190-2293	7890A/B	-
SSL Capillary Ultra Inert Inlet Low Pressure Drop Split Liner - with Glass Wool	5190-2295	7890A/B	-
PP Inlet PM kit	5188-6498	7890A/B	-
Split vent trap PM kit, single cartridge (for MMI, PTV & VI)	5188-6495	7890A/B	-
MMI Cleaning Kit	G3510-60820	7890A/B	-
PTV Septumless Head Rebuild Kit	5182-9747	7890A/B	-
PTV Septumless Head Teflon Guide	5182-9748	7890A/B	-
Ignitor (glow plug) assembly with O-ring	19231-60680	7890A/B	1
FID Collector Rebuild/Cleaning Kit	G1531-67000	7890A/B	-
Standard .011-inch FID Jet for capillary FID base	G1531-80560	7890A/B	-
High Temperature .018-inch FID Jet for capillary FID base	G1531-80620	7890A/B	-
Standard .018-inch FID Jet for packed column with packed FID base	18710-20119	7890A/B	-
Standard .011-inch FID Jet for capillary column with packed/adaptable FID base	19244-80560	7890A/B	1
High Temperature .018-inch FID Jet for capillary column with packed/adaptable FID base	19244-80620	7890A/B	-
NPD Jet, universal fit, .011-inch ID	G1534-80580	7890A/B	-
NPD Jet, universal fit, .011-inch ID Extended tip	G1534-80590	7890A/B	-
SSL Capillary Ultra Inert Inlet Gold Seal with Washer	5190-6144	7890A/B	-
SSL Capillary Ultra Inert Inlet Splitless Liner - Single taper with Glass Wool	5190-2293	7890A/B	-
**FID Collector Replacement Kit, if needed	G1531-67001	7890A/B	-

Service Engineer Comments

If there are any specific points you wish to note as part of performing the service or other items of interest for the customer, please write include them in this box.

- The Equipment can operate as normally.

Service Completion

Service request number 6007319635 Date service completed 18 Feb 2025

Agilent signature Adirek R. Customer signature _____

Total number of pages in this document 9

บริษัท แสงวิทย์ 2000 จำกัด

66/701-702 หมู่ 5 ซอยจรัญสนิทวงศ์ 13 ถนนสายบางแวก
แขวงคลองขวาง เขตภาษีเจริญ กรุงเทพฯ 10160
เลขประจำตัวผู้เสียภาษีอากร 0105543042531 สำนักงานใหญ่
โทร. (02) 861-9446-9 แฟกซ์ (02) 861-9463
E-Mail : SV2000@saengvith2000.com

**SAENGVITH 2000 CO., LTD.**

66/701-702 Moo5 Soi Charansanitwong 13, Saibangwaek Road,
Khlongkhwang, Phasicharoen, Bangkok 10160 Thailand
Tax ID : 0105543042531 Head office
Tel. (662) 861-9446-9 Fax. (662) 861-9463
E-Mail : SV2000@saengvith2000.com

TELEDYNE TEKMAR**Preventative maintenance report**

Sheet 1 of 4

Instruments Data

Repair ID:	MN-001-01-68-MN		
Instrument model:	LOTIX	serial no:	US18038002
Solid sampler:	LSS BOAT	serial no:	US18012010
Repair type:	Preventive maintenance	Date repair start:	6 February 2025
Reason for service:	Service contract	Last repair date:	8 February 2024
Electrical:	220 VAC		

Customer Data

Customer:	United Analyst and Engineering consultant	Attention:	K.Piyapat
Tel.:	02-763-2807	Fax:	02-763-2811
Address:	3 ซ.อุดมสุข 41 อ.สุขุมวิท 103 บางจาก พระโขนง กรุงเทพฯ 10260		

Software

Operating system:	Windows 10 Pro		
Teklink version:	3.0.6333.26559	Firmware version:	1.1.2169
Software communication:	USB	Communication	/ Pass Fail

Engineer:  Date: 06, 02, 25

Approved By: _____ Date: / /

เอกสารไม่ควบคุม

บริษัท แสงวิทย์ 2000 จำกัด

66/701-702 หมู่ 5 ซอยจรัญสนิทวงศ์ 13 ถนนสายบางแวก
แขวงคลองขวาง เขตภาษีเจริญ กรุงเทพฯ 10160
เลขประจำตัวผู้เสียภาษีอากร 0105543042531 สำนักงานใหญ่
โทร. (02) 861-9446-9 แฟกซ์ (02) 861-9463
E-Mail : SV2000@saengvith2000.com



SAENGVITH 2000 CO., LTD.

66/701-702 Moo5 Soi Charansanitwong 13, Saibangwaek Road,
Khlongkhwang, Phasicharoen, Bangkok 10160 Thailand
Tax ID : 0105543042531 Head office
Tel. (662) 861-9446-9 Fax. (662) 861-9463
E-Mail : SV2000@saengvith2000.com

TELEDYNE TEKMAR

Preventative maintenance report

Sheet 2 of 4

Action Data

Action	Check		
	YES	NO	Comment
Electrical: (100V/115V or 220V/240V \pm 10%)	/		219.8 V
DI water supply: Fill water reservoir and make all connections for it.	/		
Gas supply type: Free air or UHP O ₂	/		UHP
Supply gas: 500+ psi (34.5+bar) from tank to run scheduled samples	/		
Verify stage 2 pressure: 50-100 psi (3.4-6.9 bar)	/		80 psi
DI reservoir pressure (12 psi)	/		
System pressure (30-40 psi) carrier gas	/		
Carrier gas flow rate 200 cc/min \pm 10%	/		
Phosphoric acid	/		
IC sparger: Inspect & clean	/		
Mist trap: Empty then rinse the Mist trap	/		
Permeation Dryer: Inspect for damage and water accumulation	/		
Condenser loop	/		
Cooling fan	/		
Remove and inspect combustion tube and catalyst (when cooled)	/		
Clean the injection block and bottom connector	/		
Check the o-rings inside the connectors	/		
status of combustion tube	/		
Replaced catalyst	/		Check Only
Check the injection needle when cooled	/		
Sample needle / pressurization needle good?	/		
The Sample loop condition	/		
Halogen Scrubber: Check to ensure ample life for scheduled samples	/		Changed
Blue Sample filter (inline before the detector).	/		
Waste: Sufficient volume to contain the waste generate	/		

Engineer: _____

Date: 06 / 02 / 25

Approved By: _____

Date: / /

เอกสารไม่ควบคุม

บริษัท แสงวิทย์ 2000 จำกัด

66/701-702 หมู่ 5 ซอยจรัญสนิทวงศ์ 13 ถนนสายบางแวก
แขวงคลองขวาง เขตภาษีเจริญ กรุงเทพฯ 10160
เลขประจำตัวเสียภาษีอากร 0105543042531 สำนักงานใหญ่
โทร. (02) 861-9446-9 แฟกซ์. (02) 861-9463
E-Mail : SV2000@saengvith2000.com

**SAENGVITH 2000 CO., LTD.**

66/701-702 Moo5 Soi Charansanitwong 13, Saibangwaek Road,
Khlongkhwang, Phasicharoen, Bangkok 10160 Thailand
Tax ID : 0105543042531 Head office
Tel. (662) 861-9446-9 Fax. (662) 861-9463
E-Mail : SV2000@saengvith2000.com

TELEDYNE TEKMAR**Preventative maintenance report**

Sheet 3 of 4

Action Data

Diagnostics	Check		
	Pass	Fail	Comment
Injection valve	/		
Drain valve	/		
Transfer valve	/		
Needle valve	/		
IC sample valve	/		
Pressurize valve	/		
Water valve	/		
IC gas valve	/		
Sparge valve	/		
NDIR valve	/		
Carrier gas valve	/		
Mist trap valve	/		
TN Ozone valve	N/A		
Conveyor	/		
Sparger: Drain/Rinse	/		
Furnace set point	/		
Acid pump	/		
Sample needle	/		
LSS Boat injection	/		
LSS Fan/Valve	/		
Detector offset	/		

Engineer : _____

Date : 06 / 02 / 25

Approved By : _____

Date : / /

เอกสารไม่ควบคุม

บริษัท แสงวิทย์ 2000 จำกัด

66/701-702 หมู่ 5 ซอยจรัญสนิทวงศ์ 13 ถนนสายบางแวก

แขวงคลองขวาง เขตภาษีเจริญ กรุงเทพฯ 10160

เลขประจำตัวผู้เสียภาษี 0105543042531 สำนักงานใหญ่

โทร. (02) 861-9446-9 แฟกซ์. (02) 861-9463

E-Mail : SV2000@saengvith2000.com



SAENGVITH 2000 CO., LTD.

66/701-702 Moo5 Soi Charansanitwong 13, Saibangwaek Road,

Khlongkhwang, Phasicharoen, Bangkok 10160 Thailand

Tax ID : 0105543042531 Head office

Tel. (662) 861-9446-9 Fax. (662) 861-9463

E-Mail : SV2000@saengvith2000.com

TELEDYNE TEKMAR

Preventative maintenance report

Sheet 4 of 4

Action Data

Note :

ทำความสะอาด IC Sparge, Mist trap

เปลี่ยน Scrubber (Tin/Copper)

ตรวจสอบการทำงาน Heater (All), Valve (All), Conveyor, Acid pump

ตรวจสอบ Fitting & Connecting

Pressure & Flow ปกติ

เครื่องสามารถทำงานได้ตามปกติ

ควรมีอะไหล่สำรอง Scrubber (Tin/Copper), Pyrex wool, Perm dryer, Sample filter,

Ferrule 1/8 SS, 3-way Valve, Sample loop 0.5 ml, Condensate trap, Combustion tube

Catalyst, Septum 1/8 และ 1/16

Engineer: _____

Date: _____

06 02 25

Approved By: _____

Date: _____

เอกสารไม่ควบคุม



Maintenance Protocol

Atomic Fluorescence Spectrometer
mercur DUO /
mercur DUO plus

Serial-No.: K170A0153 Customer-No.: C04-006
Date: 3 February 2025 Carried out by: Mr. Srichai Fak-On

Maintenance with following Operational Qualification (OQ)
(requires a separate OQ protocol)



Company	บริษัท ยูโนเต็ด แอนนาลิสต์ แอนด์ เอ็นจิเนียริงคอนซัลแตนท์ จำกัด
User	คุณกรวิทย์
Department	ห้องปฏิบัติการ (Mercur Analysis)
Street	3 ซอยอุดมสุข 41 ถนนสุขุมวิท แขวงบางจาก เขตพระโขนง
Zip Code, City	กรุงเทพมหานคร 10260
Country	ประเทศไทย
Phone	
Fax	
E-mail	

Maintenance works basic unit

tightness visual check inside the Mercur	<input checked="" type="checkbox"/>
visual check if gold-traps are broken	<input checked="" type="checkbox"/>
visual check if spectrometer is contaminated	<input checked="" type="checkbox"/>
visual check of the fluorescence cell	<input checked="" type="checkbox"/>
visual check of the absorption cell, incl. window	<input checked="" type="checkbox"/>
reactor cleaning	<input checked="" type="checkbox"/>
check pump-hose, if necessary change it	<input checked="" type="checkbox"/>
check swivel drive (SEV)	<input checked="" type="checkbox"/>
check drying-hose, output gas-liquid-separator	<input checked="" type="checkbox"/>
test Bubble-Sensor	<input checked="" type="checkbox"/>
check gas flows	<input checked="" type="checkbox"/>
check volume flows, reagents	<input checked="" type="checkbox"/>
recording stray light values	<input checked="" type="checkbox"/>
measurement with 30 ng/l	<input checked="" type="checkbox"/>

Maintenance works Autosampler

Serial No.:

N/A

lubricate the dosing-winding (Teflon-grease-spray)	<input type="checkbox"/>
clean the dosing cylinder, if necessary exchange it	<input type="checkbox"/>
lubricate the winding system of the height drive with some drops of oil	<input type="checkbox"/>
check the toothed belt	<input type="checkbox"/>
check the position of the mechanical stopper (height: 13mm)	<input type="checkbox"/>
check the pump rate of mixing pump (<14s AS52, typ.7s/<20s AS52S, typ.10s)	<input type="checkbox"/>
check the pump rate of washing cup	<input type="checkbox"/>
check the electrical hose connections for good contact	<input type="checkbox"/>
check the connectors of the magnetic valves	<input type="checkbox"/>
check the dosing hose for buckling, if necessary exchange it	<input type="checkbox"/>

Device parameter	nominal value	actual value
visual check general tightness inside the Mercur	o.k.: <input checked="" type="checkbox"/>	changed: <input type="checkbox"/>
visual check Goldtraps	o.k.: <input checked="" type="checkbox"/>	changed: <input type="checkbox"/>
visual check spectrometer		
Fluorescence cell	o.k.: <input checked="" type="checkbox"/>	changed: <input type="checkbox"/>
Absorption cell, incl. window	o.k.: <input checked="" type="checkbox"/>	changed: <input type="checkbox"/>
lens	o.k.: <input checked="" type="checkbox"/>	changed: <input type="checkbox"/>
Swivel drive (SEV)	o.k.: <input checked="" type="checkbox"/>	changed: <input type="checkbox"/>
check pump hoses	o.k.: <input checked="" type="checkbox"/>	changed: <input type="checkbox"/>
check hoses and hose connectors	o.k.: <input checked="" type="checkbox"/>	changed: <input type="checkbox"/>
check and clean reactor	o.k.: <input checked="" type="checkbox"/>	changed: <input type="checkbox"/>
check drying hose output Gas-liquid-seperator	o.k.: <input checked="" type="checkbox"/>	changed: <input type="checkbox"/>
check bubble-sensor	o.k.: <input checked="" type="checkbox"/>	not o.k.: <input type="checkbox"/>
Check gasflow		
Argon pressure valve 4	1.2 – 1.5 bar	1.5 bar
Valve 1	10 NI/h or 0.166 NL/min	0.173 NL/min
Valve 2	50 NI/h or 0.833 NL/min	0.816 NL/min
Valve 3	5 NI/h or 0.083 NL/min	0.068 NL/min
Valve 4	10 NI/h or 0.166 NL/min	0.162 NL/min
Check liquidflow		
Acid	2.5ml/min ± 1 ml	2.5 ml/min
Red.-agent	2.5ml/min ± 1 ml	2.5 ml/min
Sample	10ml/min ± 2 ml	10 ml/min
Adventitious light - values	(V)	from file
100	0	0
200	0	0
300	0	0
350	0	0
400	0	0
450	2	1
500	4	4
550	10	9
575	14	13
600	19	18

Device parameter	nominal value	actual value
Analytical parameters Fluorescence cell		
Conditions.: max.conc.: 10µg/L PMT-voltage:451.....V		
Blank-solution		Int0.0003..
without enrichment / FBR 30 ng/L	Int > 0.0015 RSD < 3 %	Int ₁0.0028.. RSD.....1.24...%
Conditions.: max.conc.: 1.7µg/L PMT-voltage:444.....V		
Blank-solution		Int.....0.0013..
with enrichment / FBR 30 ng/L	Int > 0.008 RSD < 3 %	Int ₂0.0137.. RSD.....1.72...%
Fok.- factor (Int ₂ / Int ₁)	> 3.5	4.89
Analytical parameters Absorption cell		
Blank-solution		Ext.....0.0005..
without enrichment / FBR 100 ng/L	Ext. > 0.0012 RSD < 5 %	Ext.....0.0032.. RSD.....2.91...%
Comments		



Signature Technician

3 February 2025

Place, Date (DD/MM/YYYY)





Signature Customer

3 February 2025

Place, Date (DD/MM/YYYY)

Service Report

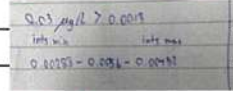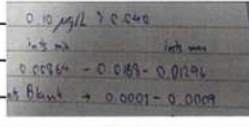Customer's address :		Customer's Ref. No.	
บริษัท ยูนิเทค แอนนาไลต์ส แอนด์ เอ็นจิเนียริงคอนซัลแตนท์ จำกัด			
3 ซอยอุดมสุข 41 ถนนสุขุมวิท แขวงบางจาก เขตพระโขนง			
กรุงเทพมหานคร 10260			
E-mail :		Phone :	Fax :
Job No. 2502070RB	User : คุณกรวิทย์	Service Engineer : ศรัชชัย พักอ่อน	Date : 3/2/2025 Page : 1/1
Instrument model : Mercury	Serial No. K170A0153	Software Version No. 4.7.10	
<input type="checkbox"/> Repair (RE) <input checked="" type="checkbox"/> Maintenance (PM) <input type="checkbox"/> Installation (IN) <input type="checkbox"/> Warranty <input type="checkbox"/> Application (AP) <input type="checkbox"/> Site Prep.(SP) <input type="checkbox"/> Visit(VI)			
Fault / Claim : Preventive maintenance Contract Year 2025 (PM 6/6)			<input type="checkbox"/> <input type="checkbox"/> <input type="checkbox"/> Error Code
Action taken : Maintenance work basic unit			
<ul style="list-style-type: none"> Check device parameters Check gas flows Check liquid flow Check adventitious light-values 		# ค่าการฟลูออเรสเซนซ์ (0.00188-0.00432) Meth: Without Enrichment. 	
Device parameter test		# ค่าการฟลูออเรสเซนซ์ (0.00864 - 0.012) Meth: Without Enrichment. 	
<ul style="list-style-type: none"> Analytical parameter fluorescence cell Analytical parameter abstraction cell 			
Action Pending / Recommendation : เครื่องมีปัญหา/ปรับ.			
<input type="checkbox"/> Spare Part <input checked="" type="checkbox"/> Instrument Configuration :			
Item No.	Name	Quantity	Unit Price
1. 407-170.240	Gas-liquid separator	1	
2. 401.580	GI clamp	1	
3. SV	PM 3 Year (2 Time/Year)	1	
4. 407.170.050	Sample inlet tube	1	
5. 407.170.057	Reactor and connector tubing 25 cm	1	
6. 407.170.069	Pre-drying membrane	1	
7. 407.170.052	Intake tubing for acid and reducing agent	2.	
Herewith the undersigned confirm the time devoted, the work performed, the perfect function of the device, and the receipt/delivery of the specified spare parts. *Traveled hours and kilometers can only be entered after the return of the service engineer.		Date / Signature of Customer	Date / Signature of Service Engineer
		สุภาวิศา	ศรัชชัย พักอ่อน
			Work completed?
			<input checked="" type="checkbox"/> Yes
			<input type="checkbox"/> No

Services are subject to the General Terms and Conditions of Analytik Jena AG, which will be sent on request.

เอกสารไม่ควบคุม

Mercur

Report file: C:\WinAAS\TMP\2025\result\WO\Pro_039
 Program version: 4.7.9.0 Printed on: 3/02/2025 13:44
 Recording started on 3/02/2025 13:32 GMT+7.0

Operator:
 Laboratory:
 Code:

Remarks:

Method parameters**Hg**

Method Without enrichment / FBR 100 ng/L PM_12-02-2024
 Created on 12/02/2024 Time 11:54
 Program ---

Parameters Mercur Technique: Hg absorption

Line	253.7 nm		
Lamp type	Hg-LP		
Integr. mode	Peak height	Integr. time	40 s
PMT	236 V		
AZ time	5 s	Peak smoothing	12/5
Delay	0 s		

Working mode	w/o enrich.	System cleaning	Acid
FBR technique	off	Wash time acid	15 s
Pump speed	4	Soaking time	20 s
Sample load time	8 s	Gas load time	10 NL/h
Reaction time	12 s		
Waiting time AZ	15 s		
Purge time1	40 s		

QC parameters

QC type	Conc. check		
QC check samp. 1	---	QC check samp. 2	---
Conc.	---	Conc.	---
Error limit	---	Error limit	---
Rep. measurement	off	Reaction	flag + continue
QC std.1 no.	1(100.00 ng/L)	QC std.2 no.	1(100.00 ng/L)
QC std.1 limit	± 50.00%	QC std.2 limit	± 0.00%
QC std. act.	flag + continue		
Expect. blank abs.	0.0100± 0.0100	Reaction	flag + continue
QC precision	off		
		Reaction	off
		QC Recal.factor	Off

Calibration settings

Calib. meth	Standard calib.	Calibr. unit	ng/L
No. standards	1	Conversion fac.	1000000
Type of standards	---	Standard prep.	Premixed
		Blank correct.	---
		Recalib. std. no.	---
Output unit	µg/L	Conversion fac.	1000
Calib. stat.	Mean	Meas. cycles	3
		Blind cycles	1
Stock sol. 1	---	Stock sol. 2	---
Stock sol. 3	---	Stock sol. 4	---
Type of cal. curve	linear	Intercept	calculated
Weighted cal.	off	Grubbs stat.	off
Check of cal. curve	no outlier test		

Sample statistics

Stat. mode	Mean	Meas. cycles	2
Confid. level	95.4 %	Blind cycles	1
Grubbs stat.	---		

Calibration standards**Hg**

No	Name	State	Pos	Conc./ ng/L	Abs	SD	RSD/%
1	Cal-Zero	(--)	##	0.00	H: 0.000544 A: 0.005800	0.000115 0.004748	21.26 81.87
2	Cal-Std1	(--)	##	100.00	H: 0.003251 A: 0.042341	0.000094 0.003312	2.921 7.824

Calibration function 1**3/02/2025 13:43 Calibration (Peak height)**


Abs=k1+k2*conc

k1=0.000544

k2=0.000027

Recal. factor:

Slope	0.00003 Abs/(ng/L)	R2-adjusted	1.0000
sc0	1.00000 ng/L	Charact. conc.	161.087 (ng/L)/1%
Lower limit	0 ng/L	Upper limit	110. ng/L
Detection limit	---	Deter. limit	---


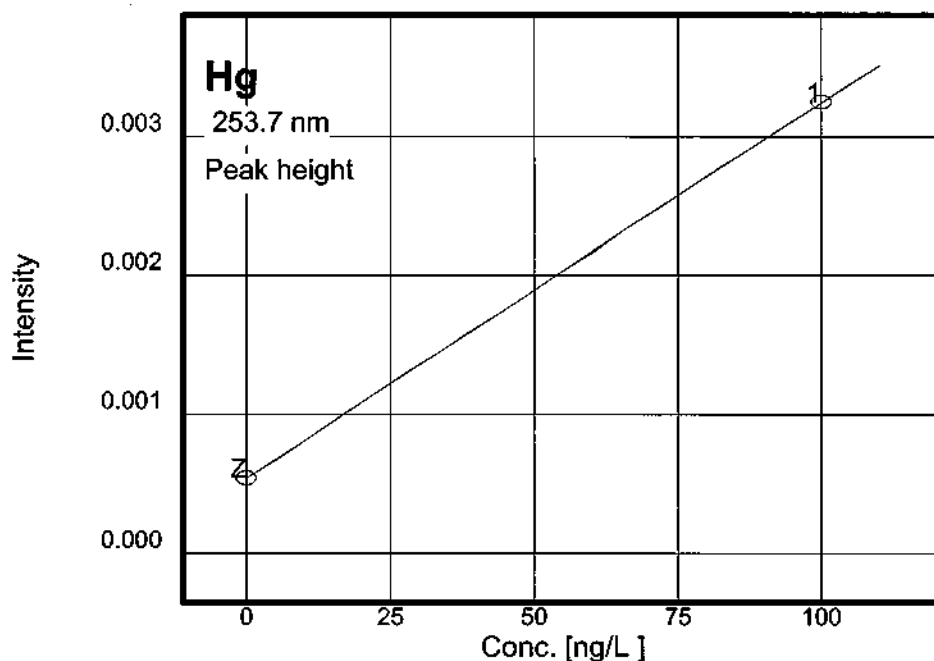


Measurements and events (sorted by time)

Hg	Without enrichment / FBR 100 ng/L PM_12-02-2024					3/02/2025	13:32
ID	Conc.	Abs	BG	SD	RSD/%	Int. type	Time
Cal-Zero		0.000564				PkH	13:35
		0.000420					13:36
		0.000649					13:37
	0ng/L	0.000544		0.00011577	21.26		13:37
Cal-Std1		0.003268				PkH	13:40
		0.003336					13:41
		0.003148					13:43
	100.ng/L	0.003251		0.000094975	2.921		13:43
Calibration	Calibration function: 01						13:43

Peak plots

Hg

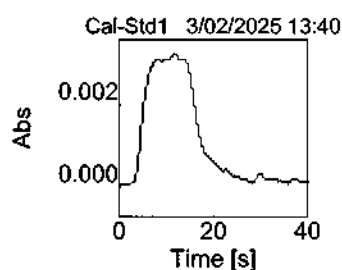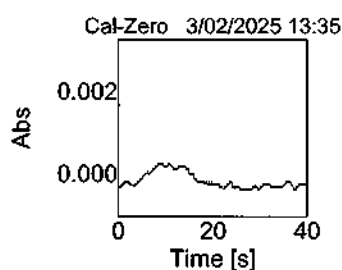

Mercur

Report file: C:\WinAAS\TMP\2025\result\WO\Pro_040
 Program version: 4.7.9.0 Printed on: 3/02/2025 17:37
 Recording started on 3/02/2025 17:26 GMT+7.0
 Operator:
 Laboratory:
 Code:

Remarks:

Method parameters**Hg**

Method Without Enrichment / FBR / 30 µg/L_PM_3-02-2025
 Created on 3/02/2025 Time 10:33
 Program ---

Parameters Mercur Technique: Hg fluorescence

Line	253.7 nm		
Lamp type	Hg-LP		
Integr. mode	Peak height	Integr. time	35 s
PMT	451 V		
AZ time	5 s	Peak smoothing	12/5
Delay	0 s		

Working mode	w/o enrich.	System cleaning	Off
FBR technique	on	Wash time acid	10 s
Pump speed	3	Soaking time	20 s
Sample load time	12 s	Gas load time	10 NL/h
Reaction time	12 s		
Waiting time AZ	5 s		
Delay	0 s		
Purge time1	30 s		
Purge time2	15 s	Gas wash time2	10 NL/h

QC parameters

QC type	Conc. check		
QC check samp. 1	---	QC check samp. 2	---
Conc.	---	Conc.	---
Error limit	---	Error limit	---
Rep. measurement	off	Reaction	flag + continue
QC std.1 no.	1(30.000 ng/L)	QC std.2 no.	3(0.100 ng/L)
QC std.1 limit	± 20.00%	QC std.2 limit	± 20.00%
QC std. act.	flag + continue		
Expect. blank abs.	0.0100± 0.0100	Reaction	flag + continue
QC precision	off		
		Reaction	off
		QC Recal.factor	Off

Calibration settings

Calib. meth	Standard calib.	Calibr. unit	ng/L
No. standards	1	Conversion fac.	1000000
Type of standards	---	Standard prep.	Premixed
		Blank correct.	---
		Recalib. std. no.	---
Output unit	µg/L	Conversion fac.	1000
Calib. stat.	Mean	Meas. cycles	3
		Blind cycles	1
Stock sol. 1	---	Stock sol. 2	---
Stock sol. 3	---	Stock sol. 4	---
Type of cal. curve	linear	Intercept	Zero
Weighted cal.	off	Grubbs stat.	off
Check of cal. curve	no outlier test		

Sample statistics

Stat. mode	Mean	Meas. cycles	3
Confid. level	95.4 %	Blind cycles	1
Grubbs stat.	off		

Calibration standards**Hg**

No	Name	State	Pos	Conc./ ng/L	Ints	SD	RSD/%
1	Cal-Zero	(-)	##	0.000	H: 0.000272 A: 0.005693	0.000004 0.000207	1.830 3.646
2	Cal-Std1	(-)	##	30.000	H: 0.002794 A: 0.03861	0.000034 0.000754	1.243 1.953

Calibration function 1

3/02/2025 17:36 Calibration (Peak height)


$$\text{Ints} = k_1 + k_2 \cdot \text{conc}$$

k1=0.000272

k2=0.000084

Recal. factor: ---

Slope	0.00008 Ints/(ng/L)	R2-adjusted	1.0000
sc0	1.00000 ng/L		
Lower limit	0 ng/L	Upper limit	33.0 ng/L
Detection limit	---	Deter. limit	---


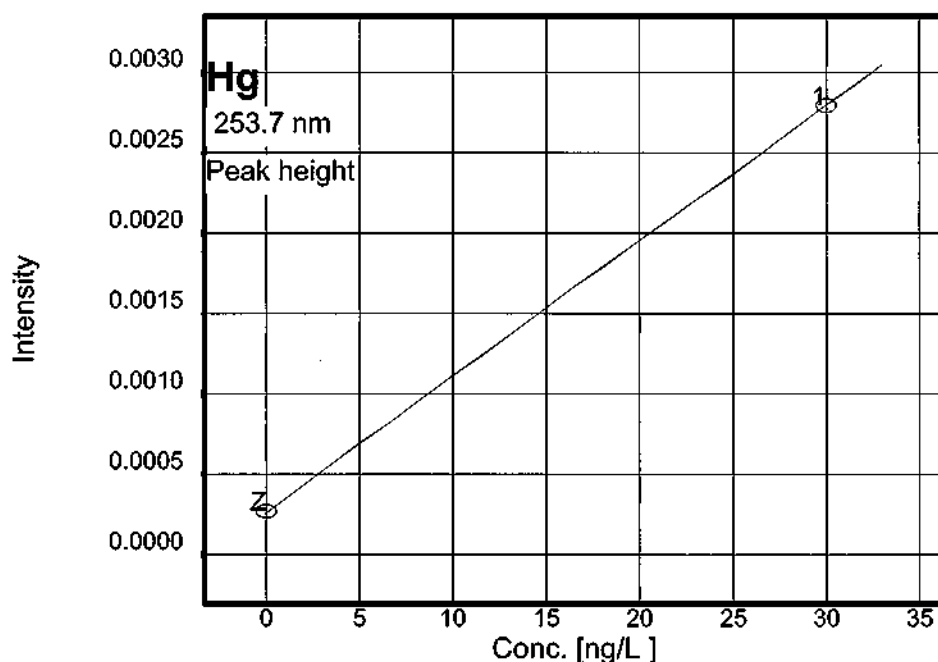


Measurements and events (sorted by time)

Hg	Without Enrichment / FBR / 30 µg/L_PM_3-02-2025					3/02/2025	17:26
ID	Conc.	Ints	BG	SD	RSD/%	Int. type	Time
Cal-Zero		0.000276				PkH	17:28
		0.000266					17:29
		0.000273					17:30
	0ng/L	0.000272		0.000004982	1.830		17:30
Cal-Std1		0.002754				PkH	17:33
		0.002812					17:34
		0.002816					17:35
	30.00ng/L	0.002794		0.000034720	1.243		17:35
Calibration	Calibration function: 01						17:36

Peak plots

Hg

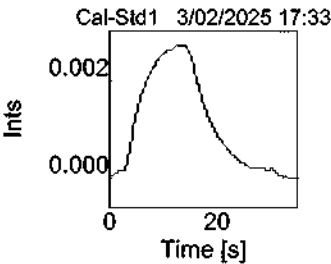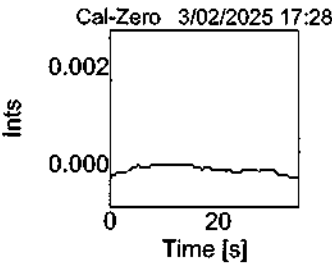

Mercur

Report file: C:\WinAAS\TMP\2025\result\WO\Pro_041
 Program version: 4.7.9.0 Printed on: 3/02/2025 17:59
 Recording started on 3/02/2025 17:46 GMT+7.0

Operator:
 Laboratory:
 Code:

Remarks:

Method parameters**Hg**

Method With Enrichment / FBR / 30 µg/L_PM_3-02-2025
 Created on 3/02/2025 Time 12:45
 Program ---

Parameters Mercur Technique: Hg fluorescence

Line	253.7 nm		
Lamp type	Hg-LP		
Integr. mode	Peak height	Integr. time	20 s
PMT	444 V		
AZ time	5 s	Peak smoothing	12/5
Delay	0 s		

Working mode	Enr. w/o reload.	System cleaning	Off
FBR technique	on	Wash time acid	10 s
Pump speed	3	Soaking time	20 s
Sample load time	10 s	Gas load time	5 NL/h
Reaction time	10 s		
Waiting time AZ	5 s		
Delay	0 s		
Purge time1	20 s		
Purge time2	15 s	Gas wash time2	5 NL/h
Purge time3	10 s	Gas wash time3	10 NL/h
Heat.time coll.1	20 s	Cool. time coll.1	25 s

QC parameters

QC type	Conc. check		
QC check samp. 1	---	QC check samp. 2	---
Conc.	---	Conc.	---
Error limit	---	Error limit	---
Rep. measurement	off	Reaction	flag + continue
QC std.1 no.	1(30.000 µg/L)	QC std.2 no.	1(30.000 µg/L)
QC std.1 limit	± 50.00%	QC std.2 limit	± 50.00%
QC std. act.	flag + continue		
Expect. blank abs.	0.0100± 0.0100	Reaction	flag + continue
QC precision	off		
		Reaction	off
		QC Recal.factor	Off

Calibration settings

Calib. meth	Standard calib.	Calibr. unit	µg/L
No. standards	1	Conversion fac.	1000
Type of standards	---	Standard prep.	Premixed
		Blank correct.	---
		Recalib. std. no.	---
Output unit	µg/L	Conversion fac.	1000
Calib. stat.	Mean	Meas. cycles	3
		Blind cycles	1
Stock sol. 1	---	Stock sol. 2	---
Stock sol. 3	---	Stock sol. 4	---
Type of cal. curve	linear	Intercept	Zero
Weighted cal.	off	Grubbs stat.	off
Check of cal. curve	no outlier test		

Sample statistics

Stat. mode	off	Meas. cycles	1
Confid. level	95.4 %	Blind cycles	1
Grubbs stat.	---		

Calibration standards**Hg**

No	Name	State	Pos	Conc./ µg/L	Ints	SD	RSD/%
1	Cal-Zero	(-)	##	0.000	H: 0.001392 A: 0.006235	0.000048 0.000289	3.475 4.635
2	Cal-Std1	(-)	##	30.000	H: 0.01371 A: 0.05663	0.000237 0.001010	1.729 1.784

Calibration function 1**3/02/2025 17:59 Calibration (Peak height)**


$$\text{Ints} = k_1 + k_2 \cdot \text{conc}$$

k1=0.001392

k2=0.000411

Recal. factor: ---



Slope	0.00041 Ints/(µg/L)	R2-adjusted	1.0000
sc0	1.00000 µg/L		
Lower limit	0 µg/L	Upper limit	33.0 µg/L
Detection limit	---	Deter. limit	---

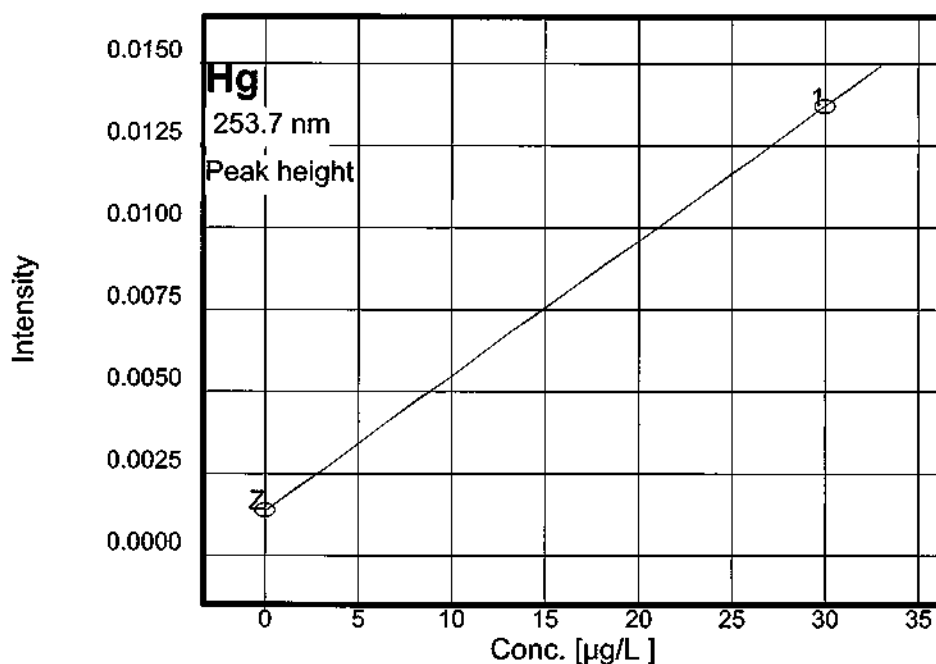**Measurements and events (sorted by time)**

Hg	With Enrichment / FBR / 30 µg/L_PM_3-02-2025					3/02/2025	17:46
ID	Conc.	Ints	BG	SD	RSD/%	Int. type	Time
Cal-Zero	0µg/L	0.001438				PkJ	17:50
		0.001397					17:51
		0.001342					17:53
		0.001392		0.000048370	3.475		17:53
Cal-Std1	30.00µg/L	0.01348				PkJ	17:56
		0.01369					17:57
		0.01395					17:59
		0.01371		0.0002370	1.729		17:59
Calibration	Calibration function: 01						17:59

Peak plots

Hg

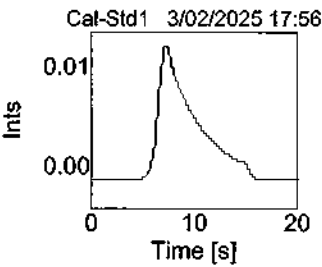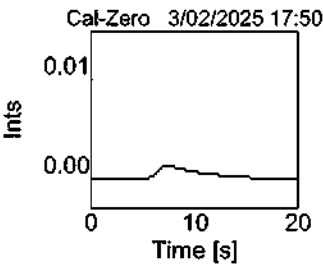

Agilent 55 240 280 Series Atomic Absorption Spectroscopy Systems

Preventive Maintenance Checklist

Agilent Preventive Maintenance provides factory recommended service for your analytical systems to assure reliable operation and the accuracy of your results.

Delivered by highly trained and certified service engineers using genuine Agilent parts and supplies, Agilent Preventive Maintenance provides everything you need to reduce unplanned downtime and keep your systems operating at their peak. This checklist will be completed at the end of the service and provided to you as a record of the installation.

Note: While non-current production AA instrument and or accessory models are not covered specifically in this document it can be used as a basic reference.

For more information about Agilent Technologies services please visit our web site using the following URL <http://www.agilent.com/en-us/services>

Introduction

Customer Information

- 1 Customers should provide all necessary operating supplies upon request of the engineer.
- 2 A customer representative should be available to the engineer while performing the preventive maintenance procedures.
- 3 Any parts, not included in the Parts Lists section of this document, are not part of the recommended Preventive Maintenance service, nor are they included in the price of this service.
- 4 If a system requires the use of extra or special procedures and/or parts for the maintenance service, then these must be ordered separately and charged as a repair, which may incur additional costs.

Important Customer Web Links

- For more information about *Agilent Technologies services*, please visit our website using the following URL: <http://www.agilent.com/en-us/products/crosslab-instrument-services/service-repair>
- To access *Agilent University*, visit <http://www.agilent.com/crosslab/university/> to learn about training options, which include online, classroom and onsite delivery. A training specialist can work directly with you to help determine your best options.
- A useful *Agilent Resource Center* web page is available, which includes short videos on maintenance, quick lists of consumables for new instruments, and other valuable information. Check out the Resource Page here: <https://www.agilent.com/en-us/agilentresources>
- Need technical support, FAQs, supplies? – visit our *Support Home page* at <http://www.agilent.com/search/support>
- Get answers. Share insights. Build connections:
Join the *Agilent Community* at <https://community.agilent.com/welcome>

Service Engineer's Responsibilities

- Contact the customer and ensure that all necessary supplies are available before the preventive maintenance visit.
- Confirm the ability of the instrument to deliver continued safe operation as established via the Agilent AA safe operation flow chart. (Refer directly to the AA 55/240/280 Preventive Maintenance Scope of Work to make this decision.)
- Only select those pages that relate to the system or module being serviced.
- Complete empty fields with the relevant information.
- Complete the relevant checkboxes in the checklist using either a "X" or tick mark "✓".
- Check "Section not applicable" check boxes to indicate services/tasks not delivered, as appropriate.
- Complete the Preventive Maintenance service in the order of the tasks listed.
- Complete the Service Review section together with the customer.
- Complete the fields for page numbers at the foot of each selected page
- Complete the total number of pages field in the Service Completion section
- Ask the customer to sign the Service Completion section including the customer's and your signature.

This information is subject to change without notice.

Instrument Maintenance

System Information

☐ Check this box if an instrument configuration report is attached instead of completing the table.

Instrument System Name and ID	240 FS AAS
Instrument System Site and Location	United Analyst and Engineering Consultant

List System Component Product Numbers	List the Serial Numbers of each Component
1. G 8432 A	M7 13160001
2.	
3.	
4.	
5.	
6.	
7.	
8.	
9.	

Preparation, Safe operation and Initial performance checks

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021

- ☐ Agilent AA safe operation flow chart inspections (to determine if the PM can be performed).

NOTE: If by following the flow chart the instrument is deemed to be unsafe for continued use you MUST NOT continue PM work. Inform the customer immediately of the Agilent recommendation that use of the instrument be discontinued.

- ☒ Discuss any specific issues with the customer before starting.
- ☐ For HF application systems, if standard sample introduction system was not installed, ask the customer to install it. *N/A*
- ☒ Review the instrument logbook for recorded problems and comments.
- ☒ Save instrument control settings before starting the procedure.
- ☒ Perform a general inspection of the system for cleanliness.
- ☒ Check for proper installation of parts, assemblies, sensors etc.
- ☒ Check system for required installation of components, settings as defined by current Service Notes
- ☒ Check for required firmware updates and verify with customers if they would like them installed.
- ☒ Use SVD to perform a Full Wavelength Scan for Cu HCL - "As found test_1"
- ☒ Perform a Basic Cu ABS test - "As found test_2"
- ☒ Print the Details page or screen captures of the test results and attach to the end of this checklist.

Preventive Maintenance Procedures

FLAME SYSTEM section

☐ *Section not applicable*

Electronic components

- ☒ Review and confirm instrument configuration data in SVD
- ☒ Confirm power supply voltages using the *SVD Power Supply diagnostic*.
- ☒ For Dual Beam instruments - Confirm RBC frequency using the *SVD RBC frequency diagnostic*.

Mechanical components

- ☒ Check the burner adjuster controls for complete and free movement. If the burner adjuster needs lubrication, use Molykote 321 or mineral-based molybdenum disulphide grease.
- ☒ Run SVD tests to exercise all motor drives over the full range of their travel:
 - ☒ Monochromator drive
 - ☒ Slit drive
 - ☒ Lamp selector
 - ☐ ABA

Optics components

- ☒ Check that external optical surfaces are clean – Clean or replace as required.
- ☒ Use SVD and perform *Mono Wavelength Correction*.
- ☒ Use SVD and perform *Slit Calibration*.
- ☒ Use SVD and perform *Grating Squareness Diagnostic*.
- ☒ Use SVD and perform *Zero Order Offset/Mono Correction*.
- ☒ Use SVD and perform *Wavelength Repeatability*.
- ☒ Physically inspect selected HC lamps (customer to supply per their choice) and measure the % Gain for each lamp. Advise customer if lamps are showing emission degradation due to age.
- ☒ Check that the signal energy of the D2 and HC lamps track properly. Advise customer if their D2 lamp is showing emission degradation due to age.

Sample Introduction and Atomization

- ☒ Inspect the burner interlock plate to ensure that the interlock pin is secure and correct for the burner type.
- ☒ Clean the burner slot with a clean white card.
- ☒ Check the uniformity of the slot width.
- ☒ Clean the burner if required.
- ☒ Change the burner o-ring.
- ☒ Clean the nebulizer, spray chamber and liquid trap.
- ☒ Change all o-rings and seals in the nebulizer, nebulizer block and spray chamber.
- ☒ Check that the pressure relief bung releases readily.
- ☒ Change o-rings on the fuel and oxidant delivery barbs
- ☒ Leave the liquid trap EMPTY and verify the flame will not ignite in this state.
- ☒ Refill liquid trap and check that overflow drains freely into the drain/waste tube.
- ☒ Check the drain/waste tube for good drainage. It should not have tight bends, kinks or loops and the lower end must be above the liquid level in the waste vessel
- ☒ Check and clean the igniter electrode

Gas handling components and safety interlocks

- ☒ Pressure test for leaks
- ☒ Leak test gasbox internal components and connections
- ☒ Check safety interlock status and operation using the *SVD interlock monitoring diagnostic*.

Analytical performance for Flame systems

- ☒ Ignite a flame.
- ☒ Check that you can adjust the nebulizer uptake rate from 4 to 6.5 mL per minute.
- ☒ Optimize the instrument ready to perform Cu sensitivity test.
- ☒ Create a manual method to perform a Basic Cu ABS test - "Final Performance Testing "
- ☒ Run a PM completed sensitivity test for a 5 ppm copper sample and record the results in the AA PM Performance test results and measurements table.

FURNACE SYSTEM section

☒ *Section not applicable*

Electronic components

- ☐ Review and confirm instrument configuration data in SVD
- ☐ Confirm power supply voltages using the *SVD Power Supply diagnostic*.

Mechanical components

- ☐ Run SVD tests to exercise all motor drives over the full range of their travel:
 - ☐ Monochromator drive
 - ☐ Slit drive
 - ☐ Lamp selector

Optics components

- ☐ Check that external optical surfaces are clean – Clean or replace as required.
- ☐ Use SVD and perform *Mono Wavelength Correction*.
- ☐ Use SVD and perform *Slit Calibration*.
- ☐ Use SVD and perform *Grating Squareness Diagnostic*.
- ☐ Use SVD and perform *Zero Order Offset/Mono Correction*.
- ☐ Use SVD and perform *Wavelength Repeatability*.
- ☐ Physically inspect selected HC lamps (customer to supply per their choice) and measure the % Gain for each lamp. Advise customer if lamps are showing emission degradation due to age.

Gas handling, water system and workhead component checks

- ☐ Inspect the GTA workhead gas hoses and connections for leaks.
- ☐ Pressure test for gas leaks
- ☐ If the cooler system is accessible (stand-alone) check for correct operation and coolant/water level – this includes any temperature and pressure settings plus filter cleaning (air flow and water).
- ☐ Inspect the GTA workhead water hoses and connections for leaks.
- ☐ Check all graphite components and replace if necessary.

Revision: 10.00, Issued: November 2021

© Agilent Technologies, Inc. 2021

- ☐ Tube
- ☐ Electrodes
- ☐ Shroud

- ☐ Check and clean the end windows on the workhead.
- ☐ Check safety interlock operation.

Analytical performance for Furnace systems

- ☐ Optimize the instrument ready to perform Cu sensitivity test.
- ☐ Run the sensitivity test for a 25 ppb copper sample and record the results in the results table.

PSD autosampler accessory for Furnace systems

☒ *Section NOT Applicable*

- ☐ Check condition of the PSD capillary – replace if necessary
- ☐ Check condition and operation of PSD syringe – ensure it does not have air locks and bubbles.
- ☐ Change PSD rinse bottle o-ring.
- ☐ Check and clean the rinse vessel.
- ☐ Check the drain tube for good drainage. It should not have tight bends, kinks or loops and the lower end must be above the liquid level in the waste vessel.
- ☐ Ensure that the waste vessel is suitable for use with the furnace system.

Sample introduction pump system (SIPS) accessory

☒ *Section NOT Applicable*

- ☐ Re-torque screws securing the hubs, presser arms and pump rotors.
- ☐ Adjust each roller so that it rotates freely.
- ☐ Wipe clean the pump rotor rollers and pump bands with a dry clean cloth.
- ☐ Ensure that the presser arms and the surfaces near the pump are free from dirt and spills.
- ☐ Remove the pump module rear cover and check for the incursion of liquids and any signs of corrosion.
- ☐ Re-torque the nuts that fasten the motor mounting plates to the chassis.
- ☐ Check clips securing the diluents holder and replace if necessary.
- ☐ Disconnect, clean T-piece, and reassemble the tubing using the following steps.

- ☐ Remove the T-piece by disconnecting the pump tubes, the pump bands and all other tubing.
- ☐ Place the T-piece in an ultrasonic bath containing strong detergent 1-5% Decon 30 or similar, for approximately 5-10 minutes.
- ☐ Wash the T-piece under a tap with a strong flow of water.
- ☐ Rinse with distilled water through all of the inlets in the reverse direction to normal sample flow.
- ☐ Reassemble.

Sample preparation system (SPS 4) accessory

☒ *Section NOT Applicable*

The Agilent SPS 4 autosampler is designed to need minimal maintenance.

The following maintenance requirements are suggested to maintain the performance of the autosampler.

- ☐ Cleaning the spill tray, rack location mat, end frames and chassis accessories with a damp soft cloth and diluted mild detergent.
- ☐ Cleaning the autosampler cover panels with domestic window cleaner.
- ☐ Checking the X- axis and Z- axis drive belts for cracks, splits, damaged teeth, excessive fraying, color changes or degradation from fumes..
- ☐ Check the X- axis, Theta- axis and Z- axis FFC cables for cracks, incorrect positioning, damaged edge or damaged connectors.

NOTE: The autosampler requires no extra lubrication throughout its lifetime.

For further details refer to the SPS 4 service manual G8410-90050.

Sample preparation system (SPS 3) accessory

☒ *Section NOT Applicable*

- ☐ Check the x-axis and z-axis timing belts – Replace if there is are any cracks, splits or color deterioration and belt tension.
- ☐ Check belt tensions - adjust if required
- ☐ Check the lubrication pad for single x-axis shaft. If pad is dry or customer has observed any vibration or erratic movements of the x-axis carriage, add 1 mL of Dow Corning 200 ® Fluid, 200 CS into the well.
- ☐ Check the auto-sampler ability to find tube positions - Calibrate if required.
- ☐ Clean the exterior surfaces of the accessory with soft lint free cloth. This cloth can be dampened with warm water or a mild detergent. Do not use organic solvents or abrasive cleaning agents.

Vapor generation accessory VGA (hydride generator)

☐ *Section NOT Applicable*

- ☐ Inspect VGA gas supply hose.
- ☐ Inspect/replace VGA pump tubing.
- ☐ Check low gas pressure interlock setting– adjust if required.
- ☐ Check precision orifice gas flow setting – adjust if required.
- ☐ Check gas regulator pressure to 46 psi (325 kPa) – adjust if required.
- ☐ Clean the exterior surfaces of the accessory with soft lint free cloth. This cloth can be dampened with warm water or a mild detergent. Do not use organic solvents or abrasive cleaning agents.

UltrAA lamp accessory (external)

☒ *Section NOT Applicable*

- ☐ Check the condition of the power cable.
- ☐ Clean the exterior surfaces of the accessory with soft lint free cloth. This cloth can be dampened with warm water or a mild detergent. Do not use organic solvents or abrasive cleaning agents.

Restore System

- ☐ If you have altered the customer's instrumentation during the course of PM, restore to the original status to allow the customer to conduct their normal activities (e.g., reload the customer's method.)

Guidance

If the PM service is performed prior to a qualification service, then use the qualification procedure as a guide for final instrument set up and checkout.

Signature Page

Service Review

- ☒ Attach available reports/printouts of all tests to this documentation.
- ☒ Record the Preventive Maintenance service activity in the customer's records/logbook.
- ☒ Update/reset instrument maintenance counters as appropriate.
- ☒ Affix the PM sticker to the system or instrument logbook based on the customer's request.
- ☒ Complete the Service Engineer Comments section if there are additional comments.
- ☒ Review this service, parts replaced, and test results obtained with the customer.
- ☒ If the instrument firmware was updated, record the details of the change in the Service Engineer's Comments box or if necessary, in the customer's IQ records.

Test Results

Test Description	Expected Test Result	Actual Test Result
Flame optics PMT Gain test		
For copper at 324.8 nm, 4 mA, 0.5 nm slit width	< 55 %	49%
Flame performance test with 5 ppm copper sample		
Air /acetylene, mixing paddle removed	Abs value > 0.5	0.5598
Air /acetylene, mixing paddle installed. 10 replicates	%RSD < 1.0	0.2%
Deuterium furnace optics PMT Gain test		
For copper at 324.8 nm, 4 mA, 0.5 nm slit width	< 55 %	-
Deuterium furnace performance test with 25 ppb copper sample (324.8 nm)		
Precision %RSD	≤ 4.0%	-
Abs value	≥ 0.15	-
Zeeman furnace analytical performance: 25 ppb copper sample (327.4 nm)		
Precision %RSD	≤ 4.0%	-
Abs value	≥ 0.10	-
MSR%	≥ 70 %	-

AA consumable and parts list table

Part Description	Part Number	Product/Model # where used	PM supplied or Consumable	Instrument-Type
Test Solution – Cu 5ppm solution	6610030100	50 55 140 240 280	PM supplied	Common
Test Solution - Blank solution	5190-7001	50 55 140 240 280	PM supplied	Common
Copper, 1000 ug/ml, 100ml	5190-8279	50 55 140 240 280	*	Common
Kit, Mk 7 O-rings, aqueous, complete set	9910093400	50 55 140 240 280	PM supplied	Flame
Organic Kit	9910093500	50 55 140 240 280	PM supplied	Flame
Wire Nebulizer Cleaning	9910024700	50 55 140 240 280	consumable	Flame
Tubing-Capillary Std Nebs	9910024800	50 55 140 240 280	consumable	Flame
Capillary Tube Hivac Neb (3) (organics only)	9910044000	50 55 140 240 280	consumable	Flame
Glass impact beads (5/pk)	9910025700	50 55 140 240 280	consumable	Flame
Teflon impact beads (5/pk): (organics only)	9910053300	50 55 140 240 280	consumable	Flame
Burner cleaning strip (100/pk)	9910053900	50 55 140 240 280	consumable	Flame
Window UV silica – round (right side)	2010082600	50 55 140 240 280	PM supplied	Common
Window UV silica – rectangular (left side)	2010082500	50 55 140 240 280	PM supplied	Common
Pad adhesive window (round)	4910012700	50 55 140 240 280	PM supplied	Common
Pad adhesive window (rectangular)	4910012800	50 55 140 240 280	PM supplied	Common
Electrode kit (1 pr) (D2)	6310003400	GTA120	PM supplied	Furnace
Shroud (D2)	6310003100	GTA120	PM supplied	Furnace
Zeeman electrode kit (1 pr)	6310003500	GTA120	PM supplied	Furnace
Zeeman shroud	6310003600	GTA120	PM supplied	Furnace
O-ring PSD rinse bottle	6910025900	PSD120	PM supplied	Furnace

* For engineers who only service AA instruments 5190-8279 can be used as a cheaper alternative for 6610030100.

Items classified as PM supplied in the above table are included in the standard PM

Those classified as consumable should be provided by the customer or charged to the customer if supplied by the Agilent service engineer.

Service Engineer Comments (optional)

If there are any specific points you wish to note as part of performing the installation or other items of interest for the customer, please write in this box.

Service Completion

Service request number 6007549143

Date service completed 30 Jan 2025

Agilent signature Kanyakorn S.

Customer signature Manida Y.

Total number of pages in this document 13

SVD Results Report



Report ID: 2 **Diagnostic Start Time:** 1/30/2025 9:14:26 AM **Diagnostic End Time:** 1/30/2025 9:46:06 AM

Customer: UAE

Service Engineer: Kanyakorn S.

Address: Soi Udomsuk 41, Sukhumvit Rd.
Bangkok

Contact Details: 026376363#1

Instrument Configuration

Configuration:

Serial Number: MY13160001	Turret Type: Automatic
Instrument Model: Varian AA140/240/280	Number Of Lamps: 4
Flame Instrument: True	Mono Type: Automatic
Furnace Instrument: True	Gasbox Type: 'Y' Gas Box
Zeeman Present: False	Auto Burner Adjuster: False
Internal Zeeman: False	Mains Frequency: 50
Internal UltraAA: False	Firmware Version: 2.11
Optics Type: Double Beam	Photomultiplier Type: Normal(900nm)
D2 BG Correction Filter: True	PWB Version: 45
Boot Block Version: 1.09	

EEPROM Data:

Instrument Run Hours: 62918.180	D2 Run Hours: 53396.500
Zero Wavelength Offset: 30.133	D2 Serial Number: not set !
Mono Correction: 0.770	D2 Install Date: 1/1/1970
Flame Hours: 32441.834	D2 Original Intensity: 1.000
	D2 Last Intensity: 475.000

Frequency:

Averaging Period: 30.0
Datapoint Count: 20

Upper Limit: 51.00	Average Frequency: 50.00	Highest Measured Frequency: 50.00
Lower Limit: 49.00		Lowest Measured Frequency: 50.00

Result: **Passed**

Power Supply:



Averaging Period: 30.0
Datapoint Count: 20

	Lower Limit (V)	Actual (V)	Upper Limit (V)	Result:
12.00 V Rail	10.80	12.12	13.20	Passed
-12.00 V Rail	-13.20	-11.90	-10.80	Passed
5.00 V Rail	4.50	5.04	5.50	Passed
310.00 V Rail	270.00	320.00	341.00	Passed

Beam Balance:

Lamp Type: Copper
Lamp Socket Used: 3

Peak Selected: 324.80

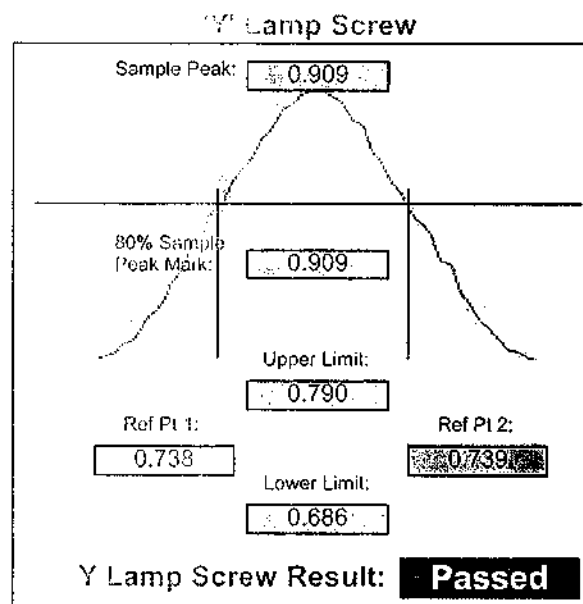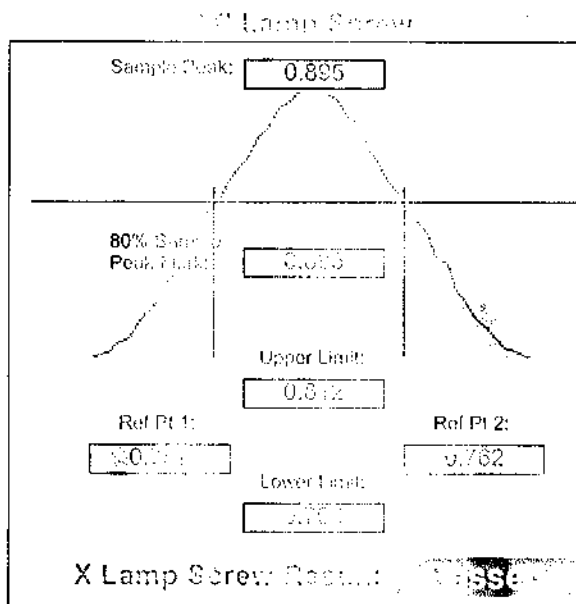Lamp Alignment: **Performed**

Grating Specifications

Lamp Element(s): Copper

Lamp Turret Position: 3

Lamp Current(mA): 4.00

3rd Width (nm): 0.5

1st Order Wavelength(nm): 324.80

Lamp Alignment: **Performed**

	Lower Limit (nm)	Central (nm)	Upper Limit (nm)	Result:
Zero Order	-0.10	0.05	0.10	Passed
First Order	324.45	324.75	325.15	Passed
Second Order	648.90	649.51	649.97	Passed

Wavelength Repeated Test:

Lamp Used: Copper	Lamp Current(mA): 4
Peak Used(nm): 324.750	Slit Width(nm): 0.2
Connected to Socket: 3	Slit Height: Normal

Lamp Alignment: 

Lower Limit(nm) 324.758	324.888 Upper Limit(nm)
-------------------------	-------------------------

(Approach from Zero Order)

(Approach from end)

Sample 1: 324.823

Sample 2: 324.823

Sample 3: 324.823

Sample 4: 324.823

Sample 5: 324.823

Sample 6: 324.819

Sample 7: 324.819


Sample 8: 324.819

Sample 9: 324.823

Sample 10: 324.819

Avg: 324.770

Standard Deviation: 0.003

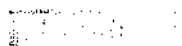Result: 

Menu:

Wavelength Drive:

Passed

Slit Drive:

Passed

Turret Drive:

Passed

Auto Burner, Monitor, Filter:

Passed

Miscellaneous:

Signal Processing Limit:

Calculate As To: New Calc Mode

	Lower Limit	Actual	Upper Limit	Result:
S0	114	1	297	Passed
S1	146	114	191	Passed
S2	271	16	332	Passed
S3	424	7	679	Passed
S4	835	11	1008	Passed
S5	1425	8	1754	Passed
S6	2408	1440	3053	Passed
S7	4747	1113	5313	Passed

Interlock:

Burner F/W: Working	Flame Detect: Working
HCO Burner F/W: Working	ECU Active: Working
Flame Shield Closed: Working	Oxidant Pressure: Working
Gas Control F/W: Working	Oxidant Changeover: Unfused
Pressure Balance Pump F/W: Working	Ignition: Working
Liquid Trap Closed: Working	

Auto Lamp Recognition:

Lamp 1: Uncoded Lamp/Not Connected	Lamp 5: Not Supported
Lamp 2: 87 - Silver/Cadmium/Lead/Zinc(UltrAA) (Ag/CLamp 6: Not Supported	
Lamp 3: 14 - Copper (Cu)	Lamp 7: Not Supported
Lamp 4: Uncoded Lamp/Not Connected	Lamp 8: Not Supported

Result: **Passed**

GTA Temperature Monitoring:

Notes:

Signatures:

UAE

Date

Kanyakorn S.
Kanyakorn S.

30 Jan 2025
Date

เอกสารไม่ควบคุม

Analyst

Date Started 1/30/2025 10:33 AM GMT: 1/30/2025 3:33 AM

Worksheet Sensitivity Test 01

Comment

Methods Cu


Computer name DESKTOP-R9UIFRS

Serial Number: MY13160001

Method: Cu (Flame)

Sample ID	Conc mg/L	%RSD	Mean Abs
CAL ZERO	0.000	38.8	0.0002
Readings			
	0.0002	0.0003	0.0001
			1/30/2025 10:51:46 AM

STANDARD 1	5.000	0.1	0.5571
Readings			
	0.5574	0.5563	0.5575
			1/30/2025 10:52:22 AM

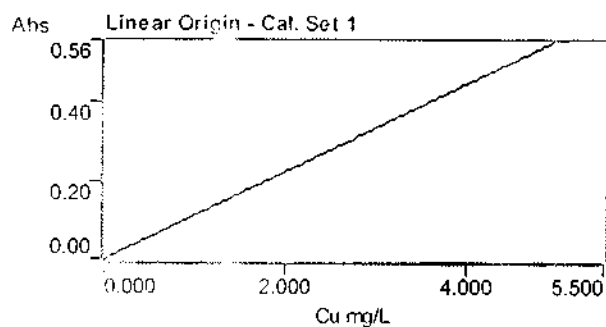

Curve Fit = Linear Origin

Characteristic Conc = 0.039 mg/L

r = 1.0000

Calculated Conc = 0.002 5 000

Residuals = -0.002 0.000

Abs = 0.11141 x C

5 ppm Cu	5.025	0.3	0.5598
Readings			
	0.5582	0.5596	0.5615
			1/30/2025 10:52:54 AM

เอกสารไม่ควบคุม

Analyst

Date Started

1/30/2025 10:33 AM GMT: 1/30/2025 3:33 AM

Worksheet

Precision Test

Comment

Methods

Cu

Computer name

DESKTOP-R9UIFRS

Serial Number:

MY13160001

Method: Cu (Flame)

Sample ID	Conc. mg/L	%RSD	Mean Abs
CAL ZERO	0.000	64.1	-0.0002

Readings


-0.0003	-0.0003	-0.0001	1/30/2025	10:46:52 AM
---------	---------	---------	-----------	-------------

STANDARD 1

5.000	0.3	0.6052
-------	-----	--------

Readings

0.6036	0.6073	0.6047	1/30/2025	10:47:24 AM
--------	--------	--------	-----------	-------------

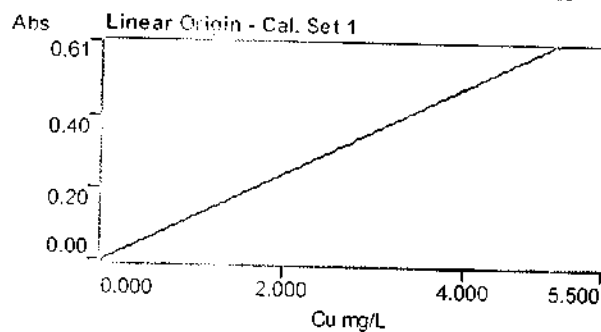

Curve Fit = Linear Origin

Characteristic Conc = 0.035 mg/L

r = 1.0000


Calculated Conc = -0.002 5.000

Residuals = 0.002 0.000


Abs = 0.12105 x C

5 ppm Cu	5.007	0.2	0.6051			
	Readings					
	0.6065	0.6052	0.6047	0.6047	0.6042	0.6079
	0.6055	0.6076	0.6064	0.6079	1/30/2025	10:48:32 AM

เอกสารไม่ควบคุม

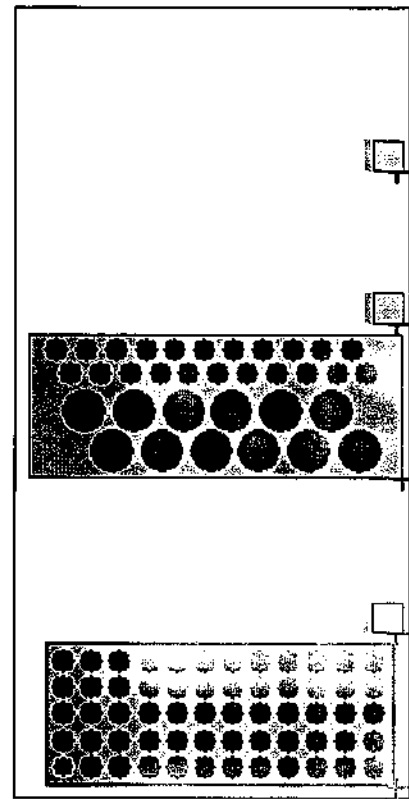

SPS 4

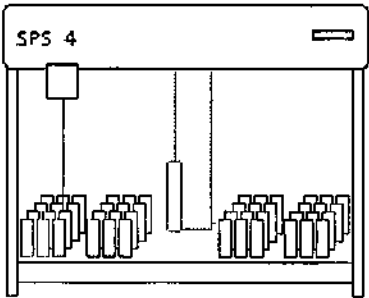

Down height (mm)

Pump speed

Key to tube colors
 Sample
 Calibration
 Calibration/QC
 Sample/QC
 Not Assigned

Sampler Offline

Goto Tube

Rack

Tube

Goto Tube

Align Probe


Rinse

Stop Rinse

Park

Optimization: Lamp

HC Lamp



0.917

Optimize Lamp

Optimize Sign

Rescale

Inst Zero

Gain 49 %

Ok

Sensitivity Check

15 mg/L gives about 0.2 Abs at 324.8 nm, A/A burner

Service Report

TO	FOR
Company: United Analyst and Engineering Consultant Co., Ltd. _ Bangkok-HQ Address: 700/2 หมู่ที่ 1 Phrakhanong District, Bangkok, 10260	Work Order Number: WO-00018067 Contact: Karnphong Boonpuang Email: karnphong.b@uaeconsultant.co.th Tel: +66 2763 2828 (7021), +66 8 6347 7390

WORK ORDER INFORMATION			
Top-Level		Order Type	Preventive Maintenance
Installed Product ID	IB-00105024	Billing Type	Paid
Product	SKALAR SAN++ Classic 2SAN59000	PO No.	SSPR2400629
Serial No.	182688	Warranty No.	
		Contract No.	

PRODUCTS SERVICED		
Installed Product Id	Serial Number	Product

PROBLEM DESCRIPTION
PM 1/1

Line Number	Engineer	Start Date And Time	End Date And Time	Billable Labor Hour	Billable Travel Hour	Travel KM
WL-00071161	Yongyuth Chanphong	02/20/2024 8:53 AM	02/20/2024 6:07 PM	9.23333		
WL-00092966	Ronnarit Dechnawarat	02/20/2024 8:53 AM	02/20/2024 6:07 PM	9.23333		
Total				18.46666	0	0

Line Number	Work Description
WL-00071161	ทำ PM เรียบร้อย
WL-00092966	ทำ PM เรียบร้อย

PARTS CONSUMED		
Part No	Part Description	Quantity

EXPENSES			
Part No	Expense Type	Description	Line Quantity


RECOMMENDED PARTS	
แนะนำอะไหล่ที่ต้องสั่งซื้อเพิ่มเติมมีดังนี้ , อะไหล่ พารามิเตอร์ Ammonia จำนวน 2 รายการ (9220, 3026) , อะไหล่ พารามิเตอร์ Phenol และ Cyanide จำนวน 6 รายการ (5454, 3028, 3031, 3034, 3036, 3150)	

REMARKS

Travel Time Disclaimer:

Please note that the travel time in this report only includes time taken to reach the installed equipment location. It does not include our engineer's return travel time.

Customer Signature:

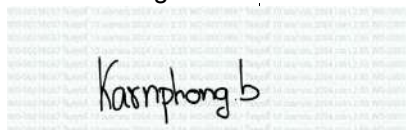

Customer Signature

Technician: Yongyuth Chanphong

Job Title: Service Manager

Email: yongyuth.yc@dksh.com

Date: 04/10/2024

Test Report

Customers	United Analyst and Engineering Consultant Co., Ltd.		
Equipment	Continuous Flow Analyzer	Manufacturer	SKALAR
Controller Mdel	SA5000	Auto Sample Model	SA1052
Controller Serial No.	182688	Auto Sample Serial No.	181729
Date of test	20-Feb-2024	Period	12 Month
Environment temperature	24.7 °C	Humidity	62.2 %RH

Results

Instrument Checked

Item	Characteristic	Before		After		Remark
1	Visual inspect	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	
2	Power supply (210 - 240 VAC)	220	VAC	220	VAC	
3	Computer	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	
4	Program	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	
5	Auto sampler	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	
6	Module holder					
	- Motor pump	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	
	- Pump tube	<input type="checkbox"/> Pass	<input checked="" type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	*
	- Air-injection	<input type="checkbox"/> Pass	<input checked="" type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	*
	- Chemistry manifolds, Switching valve, Coil, Membrane	<input type="checkbox"/> Pass	<input checked="" type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	**, ***
7	Detector					
	- Filter	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	
	- Flow cell	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	
	- Lamp	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	
8	Interface	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	
9	Rinsing valves	<input type="checkbox"/> Pass	<input type="checkbox"/> Fail	<input checked="" type="checkbox"/> N/A	<input type="checkbox"/> N/A	
10	Temperature / Reactor	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	<input type="checkbox"/> N/A
11	Flame photometer	<input type="checkbox"/> Pass	<input type="checkbox"/> Fail	<input checked="" type="checkbox"/> N/A	<input type="checkbox"/> N/A	
12	UPS / Stabilizer	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	<input checked="" type="checkbox"/> Pass	<input type="checkbox"/> Fail	<input type="checkbox"/> N/A

Warning and Error Checked

Item	Event	Before	After
13	Error list	<input checked="" type="checkbox"/> None <input type="checkbox"/> Appear : _____	<input checked="" type="checkbox"/> None <input type="checkbox"/> Appear : _____

Check with Standard

Item	Characteristic	Before	After	Remark
14	Base Line Test	<input type="checkbox"/> Pass <input type="checkbox"/> Fail <input checked="" type="checkbox"/> N/A	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail <input type="checkbox"/> N/A	
15	Detector Signal Test	<input type="checkbox"/> Pass <input type="checkbox"/> Fail <input checked="" type="checkbox"/> N/A	<input checked="" type="checkbox"/> Pass <input type="checkbox"/> Fail <input type="checkbox"/> N/A	

Summary of checked

- ☒ The instrument can work normally and efficiently. (เครื่องมือวัดสามารถทำงานได้ปกติและมีประสิทธิภาพ)
- ☐ The instrument can work but it's requiring to maintenance. (เครื่องมือวัดสามารถทำงานได้แต่ต้องบำรุงรักษา)
- ☐ The instrument could not work it's requiring to repair. (เครื่องมือวัดไม่สามารถทำงานได้แต่ต้องการซ่อมบำรุง)

Remark :

* Pump tube และ Air tube เริ่มเสื่อมสภาพ ได้เปลี่ยนอะไหล่ทั้งหมดแล้วตามระยะการใช้งาน
** อะไหล่ 9220 (Manifold T ,StSt needle) ขำรดไม่สามารถซ่อมได้ และได้เปลี่ยนมาใช้ 5216 ทดแทนชั่วคราว
*** อะไหล่ 5454 (Nipple polyethylene N5) เสีย ได้ทำการเปลี่ยนใหม่แล้ว
หมายเหตุ แนะนำอะไหล่ที่ควรสั่งเพิ่มเติมดังนี้
1. อะไหล่ พารามิเตอร์ Ammonia จำนวน 2 รายการ (9220, 3026)
2. อะไหล่ พารามิเตอร์ Phenol และ Cyanide จำนวน 6 รายการ (5454, 3028, 3031, 3034, 3036, 3150)

Standard Equipment Used

Equipment	Equipment I.D.	
Digital multi meter	S/N 57600592	Due date : 8-Jul-2024
Thermo hygrometer	S/N 39520444	Due date : 27-Dec-2024

Test By :

(Mr. Yongyuth Chanphong)

Approved by :


(Mr. Eknamong Wankliang)


Position :


Supervisor, Technical Service

Position :


Manager, Technical Services

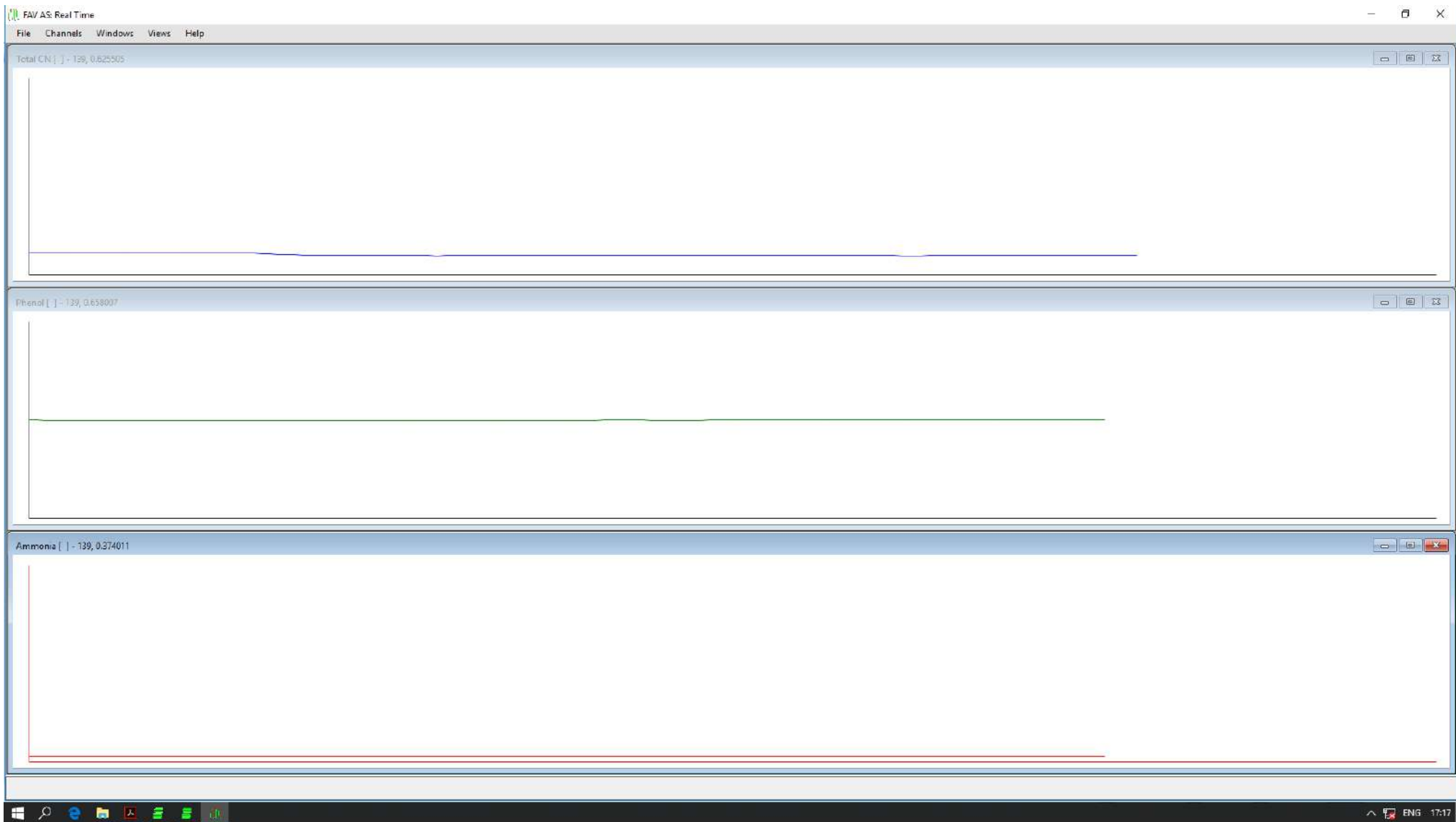

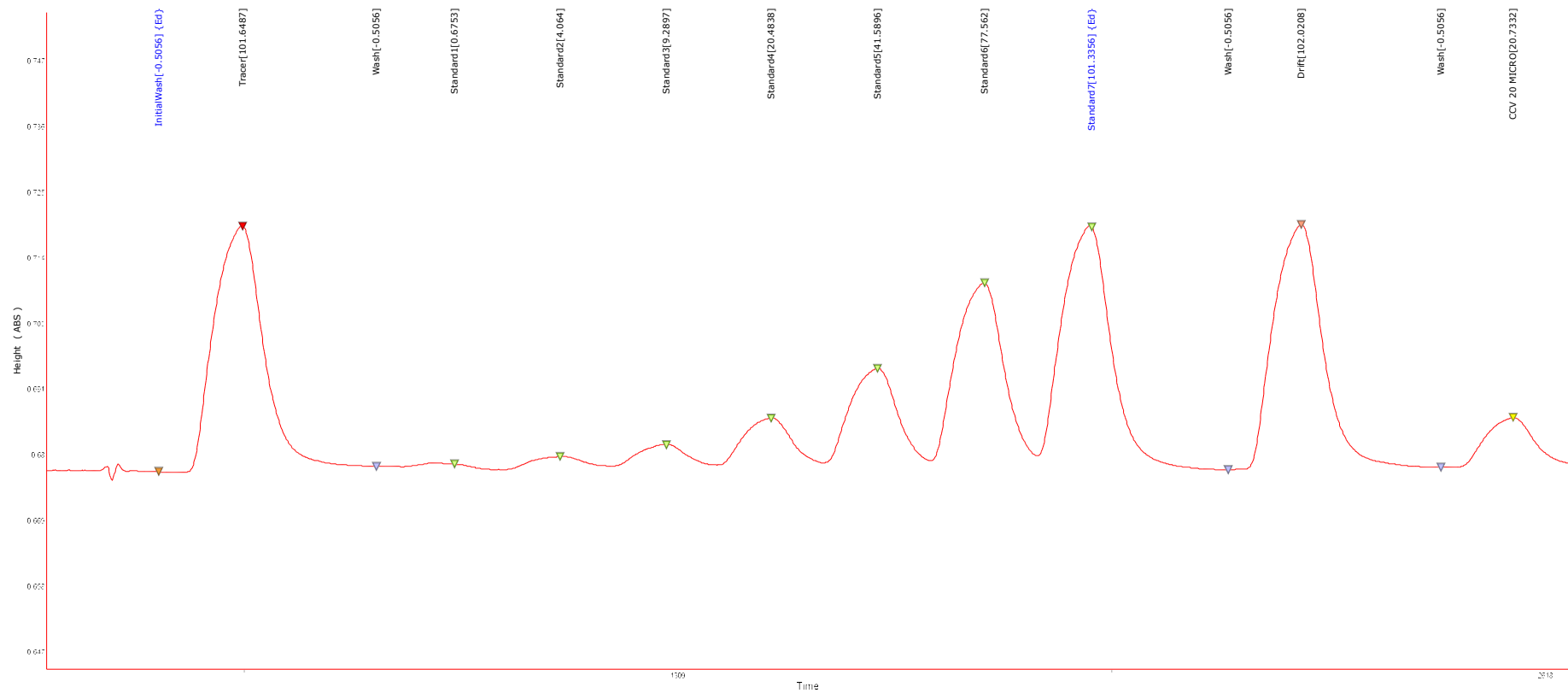

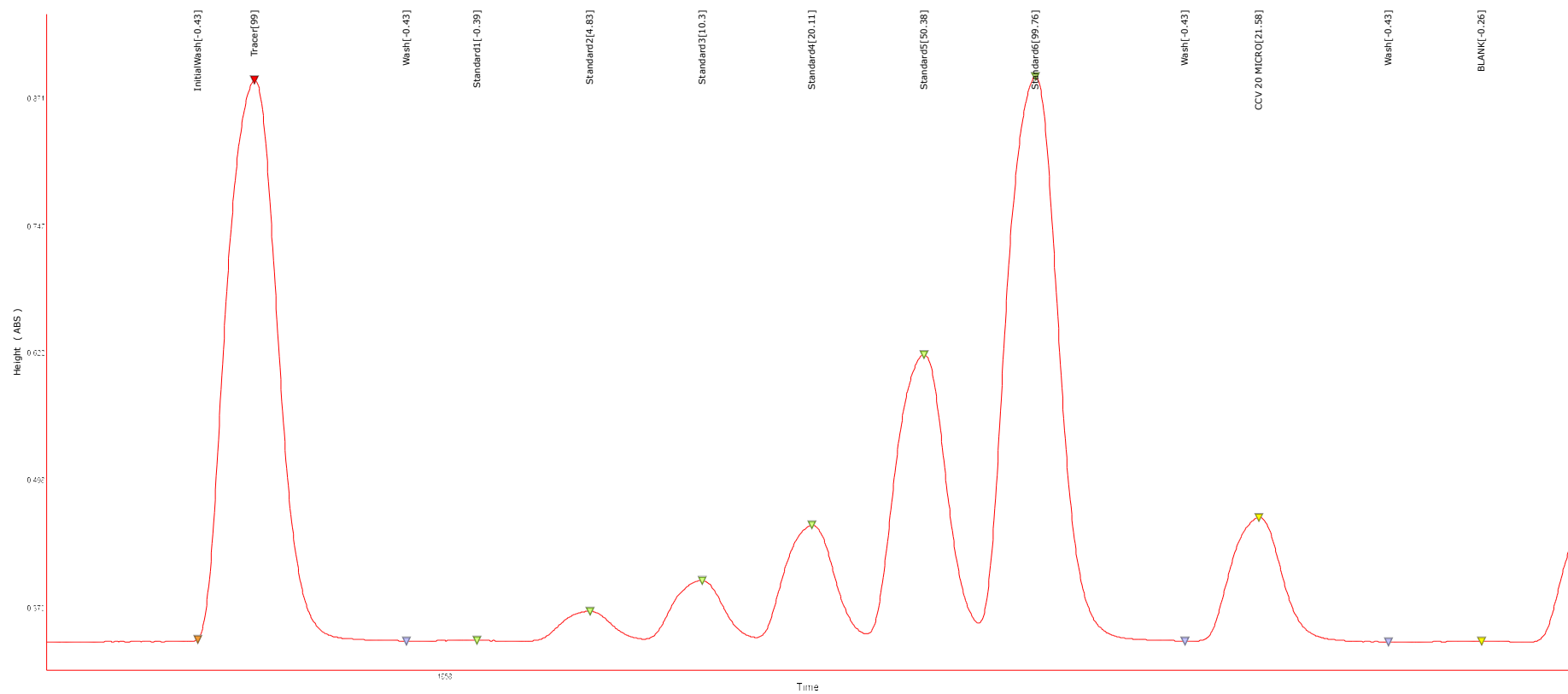

Total CN
FlowAccess V3 : Date : 19/04/2024 16:07
Run Name : CN20240222A1 - Run Database Reference Name : CN20240222A1S
User Name : Administrator Operator Name : Administrator

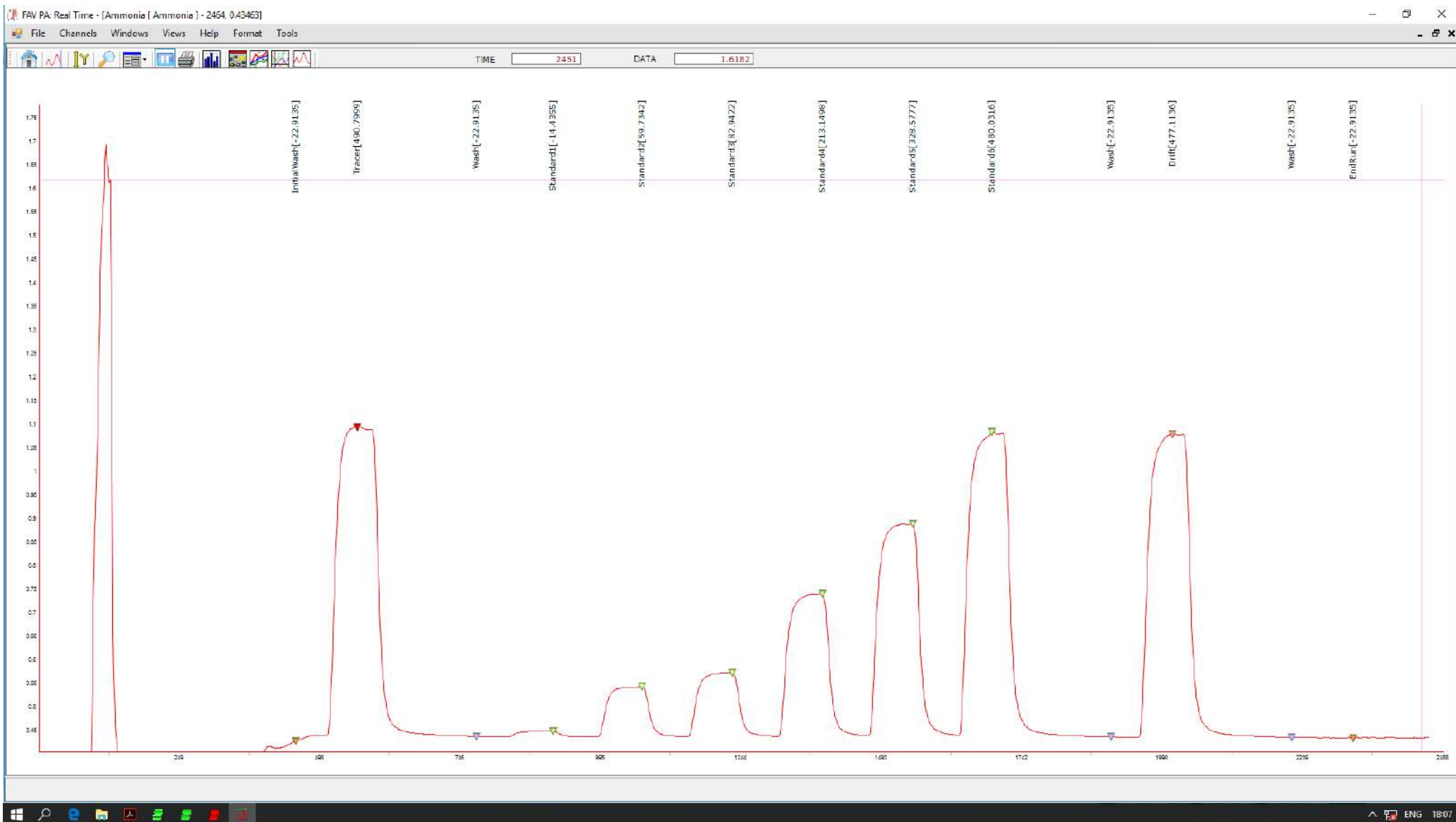

CERTIFICATE OF CALIBRATION

Certificate No. : SP24-001

Page 1 of 5

Customer : United Analyst and Engineering Consultant Co.,Ltd. (Head Office)

Address : 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260

Location of calibration : Laboratory 213

Equipment : UV-Vis Spectrophotometer

Manufacturer : Hitachi

Model : U-2900

Serial No. : 21E22-009

ID No. : UAE.WAT.051/2564

Received Date : 4 January 2024

Calibration Date : 4 January 2024

Issue Date : 5 January 2024


Condition Instrument : Good

Calibrated by :


(Mr.Tanawut Rittidach)

Technical Manager

Approved by :

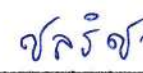
(Ms. Chonthicha Sangngern)

Quality Manager

The calibration result is applied only to the above calibrated item and was found accurate as shown on date and place of calibration only.

The measurement capability of the laboratory and its traceability to recognized national standards and to the unit of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the DQE Services Co., Ltd.

เอกสารไม่ควบคุม

REPORT OF CALIBRATION

Certificate No. : SP24-001

Page 2 of 5

Environment Condition : Ambient Temperature 25 ± 5 °CRelative humidity 55 ± 20 %RH**Calibration method :** In-house method CP-01 Based on ASTM E275-08**Certified Reference Materials :**

Material	Serial No.	Certificate No.	Due date
Absorbance Standard set	25760	115663	25 October 2025
Absorbance Standard set	25757	115638	25 October 2025
Wavelength Standard set	25806	115657	25 October 2025
Wavelength Standard set	25758	115665	25 October 2025

Traceability : This certification is traceable to the International System of Unit maintained at National -
Institute of Standards and Technology (NIST) through Starna Scientific Limited

Spectral Band Width of UUC : 1.5 nm.**Scan Speed of UUC :** 200 nm/min**Scan Interval of UUC :** 0.1 nm.**Resolution of UUC :** Photometric 0.001 Abs.

Wavelength 0.1 nm.

REPORT OF CALIBRATION

Certificate No. : SP24-001

Page 3 of 5

Calibration Results : Without adjustment

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor <i>k</i>
420	0.0000	0.000	0.0000	0.0028	2.00
	0.5780	0.575	0.0030	0.0031	2.00
	1.0484	1.045	0.0034	0.0029	2.00
	2.1876	2.192	-0.0044	0.0080	2.00
440	0.0000	0.000	0.0000	0.0028	2.00
	0.5595	0.558	0.0015	0.0034	2.00
	1.0239	1.023	0.0009	0.0035	2.00
	2.1230	2.125	-0.0020	0.0079	2.00
465	0.0000	0.000	0.0000	0.0028	2.00
	0.5230	0.520	0.0030	0.0030	2.00
	0.9633	0.961	0.0023	0.0029	2.00
	1.9753	1.975	0.0003	0.0070	2.00
546.1	0.0000	0.000	0.0000	0.0028	2.00
	0.5181	0.516	0.0021	0.0031	2.00
	1.0002	0.997	0.0032	0.0033	2.00
	1.9973	1.993	0.0043	0.0084	2.00
590	0.0000	0.000	0.0000	0.0028	2.00
	0.5517	0.550	0.0017	0.0030	2.00
	1.0803	1.079	0.0013	0.0030	2.00
	2.0373	2.032	0.0053	0.0080	2.00
635	0.0000	0.000	0.0000	0.0028	2.00
	0.5591	0.558	0.0011	0.0031	2.00
	1.0518	1.050	0.0018	0.0030	2.00
	1.9274	1.923	0.0044	0.0079	2.00

เอกสารไม่ควบคุม

REPORT OF CALIBRATION

Certificate No. : SP24-001

Page 4 of 5

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor <i>k</i>
235	0.0000	0.000	0.0000	0.0050	2.00
	0.7469	0.743	0.0039	0.0057	2.00
257	0.0000	0.000	0.0000	0.0050	2.00
	0.8674	0.862	0.0054	0.0059	2.00
313	0.0000	0.000	0.0000	0.0050	2.00
	0.2919	0.289	0.0029	0.0051	2.00
350	0.0000	0.000	0.0000	0.0050	2.00
	0.6430	0.641	0.0020	0.0055	2.00

REPORT OF CALIBRATION

Certificate No. : SP24-001

Page 5 of 5

Wavelength Accuracy :

CRMs Values (nm.)	UUC Reading (nm.)	Correction (nm.)	Uncertainty (nm.)	Coverage factor <i>k</i>
241.72	241.2	0.52	0.18	2.00
279.45	279.0	0.45	0.18	2.00
287.81	287.4	0.41	0.18	2.00
334.06	333.8	0.26	0.18	2.00
360.93	360.6	0.33	0.18	2.00
418.59	418.4	0.19	0.18	2.00
445.94	445.8	0.14	0.18	2.00
453.66	453.4	0.26	0.18	2.00
460.02	459.8	0.22	0.18	2.00
536.59	536.4	0.19	0.18	2.00
637.98	638.0	-0.02	0.18	2.00
431.38	431.2	0.18	0.18	2.00
472.50	472.5	0.00	0.18	2.00
513.47	513.4	0.07	0.18	2.00
528.88	528.9	-0.02	0.18	2.00
573.17	573.4	-0.23	0.18	2.00
585.35	585.2	0.15	0.20	2.00
684.40	684.4	0.00	0.18	2.00
740.72	741.0	-0.28	0.20	2.00
748.55	748.8	-0.25	0.18	2.00
807.03	807.1	-0.07	0.18	2.00
879.28	879.5	-0.22	0.18	2.00

Remark : - UUC = Unit Under Calibration

- N/A = Not Available

- The result expanded uncertainty of measurement U is stated as the standard uncertainty of measurement multiplied by the coverage factor k ,

which for a normal distribution corresponds to a coverage probability of approximately 95%

- * Indicates non TISI accredited


- End of Certificate -

เอกสารไม่ควบคุม

CERTIFICATE OF CALIBRATION

Certificate No. : SP24-028

Page 1 of 5

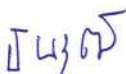Customer : United Analyst and Engineering Consultant Co.,Ltd. (Head Office)**Address :** 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260**Location of calibration :** Laboratory 315**Equipment :** UV-Vis Spectrophotometer**Manufacturer :** HITACHI**Model :** U-5100**Serial No. :** 23A4-008**ID No. :** UAE.WAS.010/2567**Received Date :** 10 September 2024**Calibration Date :** 10 September 2024**Issue Date :** 13 September 2024**Condition Instrument :** Good**Calibrated by :**

(Mr.Tanawut Rittidach)

Technical Manager

Approved by :

(Ms. Chonthicha Sangngern)

Quality Manager

The calibration result is applied only to the above calibrated item and was found accurate as shown on date and place of calibration only.

The measurement capability of the laboratory and its traceability to recognized national standards and to the unit of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the DQE Services Co., Ltd.

เอกสารไม่ควบคุม

REPORT OF CALIBRATION

Certificate No. : SP24-028

Page 2 of 5

Environment Condition : Ambient Temperature 25 ± 5 °CRelative humidity 55 ± 20 %RH**Calibration method :** In-house method CP-01 Based on ASTM E275-08**Certified Reference Materials :**

Material	Serial No.	Certificate No.	Due date
Absorbance Standard set	25760	115663	25 October 2025
Absorbance Standard set	25757	115638	25 October 2025
Wavelength Standard set	25806	115657	25 October 2025
Wavelength Standard set	25758	115665	25 October 2025

Traceability : This certification is traceable to the International System of Unit maintained at National -
Institute of Standards and Technology (NIST) through Sarna Scientific Limited

Spectral Band Width of UUC : 5.0 nm.**Scan Speed of UUC :** 40**Scan Interval of UUC :** 0.1 nm.**Resolution of UUC :** Photometric 0.001 Abs.

Wavelength 0.1 nm.

REPORT OF CALIBRATION

Certificate No. : SP24-028

Page 3 of 5

Calibration Results : Without adjustment

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor <i>k</i>
420	0.0000	0.000	0.0000	0.0028	2.00
	0.5780	0.575	0.0030	0.0031	2.00
	1.0484	1.044	0.0044	0.0029	2.00
	2.1876	2.190	-0.0024	0.0075	2.00
440	0.0000	0.000	0.0000	0.0028	2.00
	0.5595	0.557	0.0025	0.0034	2.00
	1.0239	1.021	0.0029	0.0035	2.00
	2.1230	2.121	0.0020	0.0079	2.00
465	0.0000	0.000	0.0000	0.0028	2.00
	0.5230	0.519	0.0040	0.0029	2.00
	0.9633	0.961	0.0023	0.0028	2.00
	1.9753	1.975	0.0003	0.0070	2.00
546.1	0.0000	0.000	0.0000	0.0028	2.00
	0.5181	0.515	0.0031	0.0031	2.00
	1.0002	0.997	0.0032	0.0033	2.00
	1.9973	1.996	0.0013	0.0085	2.00
590	0.0000	0.000	0.0000	0.0028	2.00
	0.5517	0.549	0.0027	0.0030	2.00
	1.0803	1.078	0.0023	0.0029	2.00
	2.0373	2.031	0.0063	0.0081	2.00
635	0.0000	0.000	0.0000	0.0028	2.00
	0.5591	0.557	0.0021	0.0031	2.00
	1.0518	1.049	0.0028	0.0029	2.00
	1.9274	1.923	0.0044	0.0080	2.00

REPORT OF CALIBRATION

Certificate No. : SP24-028

Page 4 of 5

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor <i>k</i>
235	0.0000	0.000	0.0000	0.0050	2.00
	0.7469	0.743	0.0039	0.0056	2.00
257	0.0000	0.000	0.0000	0.0050	2.00
	0.8674	0.862	0.0054	0.0059	2.00
313	0.0000	0.000	0.0000	0.0050	2.00
	0.2919	0.291	0.0009	0.0051	2.00
350	0.0000	0.000	0.0000	0.0050	2.00
	0.6430	0.639	0.0040	0.0055	2.00

REPORT OF CALIBRATION

Certificate No. : SP24-028

Page 5 of 5

Wavelength Accuracy :

CRMs Values (nm.)	UUC Reading (nm.)	Correction (nm.)	Uncertainty (nm.)	Coverage factor <i>k</i>
241.00	240.4	0.60	0.18	2.00
279.30	278.7	0.60	0.18	2.00
288.90	288.5	0.40	0.18	2.00
334.50	334.2	0.30	0.18	2.00
361.40	361.1	0.30	0.18	2.00
418.40	418.0	0.40	0.18	2.00
447.20	446.7	0.50	0.18	2.00
459.30	459.6	-0.30	0.18	2.00
537.00	536.6	0.40	0.18	2.00
638.00	637.4	0.60	0.18	2.00
441.29	440.8	0.49	0.18	2.00
479.88	479.6	0.28	0.18	2.00
513.75	513.5	0.25	0.18	2.00
528.59	528.6	-0.01	0.18	2.00
575.10	574.9	0.20	0.18	2.00
585.56	585.3	0.26	0.20	2.00
684.70	684.1	0.60	0.18	2.00
740.51	740.0	0.51	0.20	2.00
747.61	747.2	0.41	0.18	2.00
807.04	806.3	0.74	0.18	2.00
879.68	878.9	0.78	0.18	2.00

Remark : - UUC = Unit Under Calibration

- N/A = Not Available

- The result expanded uncertainty of measurement U is stated as the standard uncertainty of measurement multiplied by the coverage factor k ,


which for a normal distribution corresponds to a coverage probability of approximately 95%

- End of Certificate -

CERTIFICATE OF CALIBRATION


Certificate No. : SP25-001

Page 1 of 5

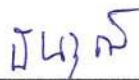Customer : United Analyst and Engineering Consultant Co.,Ltd. (Head Office)**Address :** 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260**Location of calibration :** Laboratory 213**Equipment :** UV-Vis Spectrophotometer**Manufacturer :** Hitachi**Model :** U-2900**Serial No. :** 21E22-009**ID No. :** UAE.WAT.051/2564**Received Date :** 3 January 2025**Calibration Date :** 3 January 2025**Issue Date :** 8 January 2025**Condition Instrument :** Good**Calibrated by :**

(Mr.Tanawut Rittidach)

Technical Manager

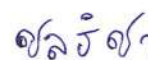Approved by :

(Ms. Chonthicha Sangngern)

Quality Manager

The calibration result is applied only to the above calibrated item and was found accurate as shown on date and place of calibration only.

The measurement capability of the laboratory and its traceability to recognized national standards and to the unit of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the DQE Services Co., Ltd.

เอกสารไม่ควบคุม

REPORT OF CALIBRATION

Certificate No. : SP25-001

Page 2 of 5

Environment Condition : Ambient Temperature 25 ± 5 °CRelative humidity 55 ± 20 %RH**Calibration method :** In-house method CP-01 Based on ASTM E275-08**Certified Reference Materials :**

Material	Serial No.	Certificate No.	Due date
Absorbance Standard set	25760	115663	25 October 2025
Absorbance Standard set	25757	115638	25 October 2025
Wavelength Standard set	25806	115657	25 October 2025
Wavelength Standard set	25758	115665	25 October 2025

Traceability : This certification is traceable to the International System of Unit maintained at National -

Institute of Standards and Technology (NIST) through Starna Scientific Limited

Spectral Band Width of UUC : 1.5 nm.**Scan Speed of UUC :** 200 nm/min**Scan Interval of UUC :** 0.1 nm.**Resolution of UUC :** Photometric 0.001 Abs.

Wavelength 0.1 nm.

REPORT OF CALIBRATION

Certificate No. : SP25-001

Page 3 of 5

Calibration Results : Without adjustment

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor <i>k</i>
420	0.0000	0.000	0.0000	0.0028	2.00
	0.5780	0.578	0.0000	0.0031	2.00
	1.0484	1.045	0.0034	0.0029	2.00
	2.1876	2.192	-0.0044	0.0075	2.00
440	0.0000	0.000	0.0000	0.0028	2.00
	0.5595	0.560	-0.0005	0.0034	2.00
	1.0239	1.023	0.0009	0.0035	2.00
	2.1230	2.125	-0.0020	0.0079	2.00
465	0.0000	0.000	0.0000	0.0028	2.00
	0.5230	0.521	0.0020	0.0030	2.00
	0.9633	0.961	0.0023	0.0029	2.00
	1.9753	1.977	-0.0017	0.0070	2.00
546.1	0.0000	0.000	0.0000	0.0028	2.00
	0.5181	0.518	0.0001	0.0031	2.00
	1.0002	0.998	0.0022	0.0033	2.00
	1.9973	1.993	0.0043	0.0084	2.00
590	0.0000	0.000	0.0000	0.0028	2.00
	0.5517	0.552	-0.0003	0.0030	2.00
	1.0803	1.079	0.0013	0.0030	2.00
	2.0373	2.032	0.0053	0.0079	2.00
635	0.0000	0.000	0.0000	0.0028	2.00
	0.5591	0.559	0.0001	0.0031	2.00
	1.0518	1.050	0.0018	0.0030	2.00
	1.9274	1.923	0.0044	0.0079	2.00

เอกสารไม่ควบคุม

REPORT OF CALIBRATION

Certificate No. : SP25-001

Page 4 of 5

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor <i>k</i>
235	0.0000	0.000	0.0000	0.0050	2.00
	0.7469	0.744	0.0029	0.0057	2.00
257	0.0000	0.000	0.0000	0.0050	2.00
	0.8674	0.863	0.0044	0.0059	2.00
313	0.0000	0.000	0.0000	0.0050	2.00
	0.2919	0.290	0.0019	0.0051	2.00
350	0.0000	0.000	0.0000	0.0050	2.00
	0.6430	0.640	0.0030	0.0055	2.00

REPORT OF CALIBRATION

Certificate No. : SP25-001

Page 5 of 5

Wavelength Accuracy :

CRMs Values (nm.)	UUC Reading (nm.)	Correction (nm.)	Uncertainty (nm.)	Coverage factor <i>k</i>
241.72	241.1	0.62	0.18	2.00
279.45	279.0	0.45	0.18	2.00
287.81	287.3	0.51	0.18	2.00
334.06	333.8	0.26	0.18	2.00
360.93	360.6	0.33	0.18	2.00
418.59	418.2	0.39	0.18	2.00
445.94	445.5	0.44	0.18	2.00
453.66	453.4	0.26	0.18	2.00
460.02	459.8	0.22	0.18	2.00
536.59	536.6	-0.01	0.18	2.00
637.98	637.7	0.28	0.18	2.00
431.38	431.1	0.28	0.18	2.00
472.50	472.3	0.20	0.18	2.00
513.47	513.4	0.07	0.18	2.00
528.88	528.9	-0.02	0.18	2.00
573.17	573.3	-0.13	0.18	2.00
585.35	585.1	0.25	0.20	2.00
684.40	684.5	-0.10	0.18	2.00
740.72	741.0	-0.28	0.20	2.00
748.55	748.8	-0.25	0.18	2.00
807.03	807.3	-0.27	0.18	2.00
879.28	879.6	-0.32	0.18	2.00

Remark : - UUC = Unit Under Calibration

- N/A = Not Available

- The result expanded uncertainty of measurement U is stated as the standard uncertainty of measurement multiplied by the coverage factor k ,

which for a normal distribution corresponds to a coverage probability of approximately 95%

- End of Certificate -

เอกสารไม่ควบคุม



THAI UNIQUE OPEN LAB SERVICE

OPERATIONAL QUALIFICATION REPORT (OQ)

Equipment Operational Qualification Report

Report No.	SV2407/21898
Equipment	GC-MS
System Model	SQ
System ID	GQS1203F021
Equipment Details	Bruker
Test Protocol	Scion OQ Protocol
Protocol Rev.	A
Date	16-Jul-24
Report Type	Report
Org. Name	United Analyst and Engineering Consultant Co.,Ltd
Org. Location	3 Soi Udomsuk 41 Sukhumvit Rd. Bangchak Phrakhanong Bagkok Thailand 10260

เอกสารไม่ควบคุม



บริษัท ไทยยูนิค จำกัด

THAI UNIQUE CO., LTD.

80-82 ถนนประชาธิปไตย แขวงบางขุนพรหม เขตพระนคร กรุงเทพฯ 10200

80-82 Prachathipatai Rd., Bangkhunphrom, Pranakorn, Bangkok 10200

Tel. 0-2629-0191-6, 0-2280-1787, Fax. 0-2280-1788, E-mail : thawatt@thaiunique.com, Website : www.thaiunique.com

CERTIFICATE OF CALIBRATION

GAS CHROMATOGRAPH MASS SPECTROMETER

Certificate No.: SV2407/21898

Customer: United Analyst and Engineering Consultant Co., Ltd.

Address: 3 Soi Udomsuk 41 Sukhumvit Rd. Bangchak
Phrakhanong Bagkok Thailand 10260

Instruments Model:	MS Scion-SQ	S/N	GQS1203F021
	GC 451-GC	S/N	BR1203M099
	AUTO SAMPLER CP8400	S/N	BR1203M331

Standard Reference Number: 393065201

Procedure Document Number: 394207000

System Test

PM perform and Diagnostic Test	<input checked="" type="checkbox"/> PASS	<input type="checkbox"/> FAIL
Air Water Check Test	<input checked="" type="checkbox"/> PASS	<input type="checkbox"/> FAIL
Tune Test EI	<input checked="" type="checkbox"/> PASS	<input type="checkbox"/> FAIL
Signal to Noise Test (EI) SCAN	<input checked="" type="checkbox"/> PASS	<input type="checkbox"/> FAIL
Injection EI Area Precision Test	<input checked="" type="checkbox"/> PASS	<input type="checkbox"/> FAIL
Injection EI RT Precision Test	<input checked="" type="checkbox"/> PASS	<input type="checkbox"/> FAIL
User Demonstration	<input checked="" type="checkbox"/> PASS	<input type="checkbox"/> FAIL

Engineer


Somchai Pohtongkam

Date

16 July 2024



Thai Unique Co., Ltd.

Service Division

เอกสารไม่ควบคุม

SCION™

Operational Qualification Protocol

For SCION Instrument

Name and Model:

Serial Number:

System ID Number:

Publication no. 394207000

Revision A

November 2011

Contact

Scion Customer Service and Support uses a Customer Relationship Management (CRM) system. The interaction with this system offers the Customer immediate benefits including the contact center or help desk.

Scion worldwide service & support offices can be found from Scion website:



www.scion.com/support.html

Scion is the owner of copyright on this document and any associated software. Under law, the written permission of Scion must be obtained before the documentation or the software is copied, reproduced, translated or converted to electronic or other machine-readable form, in whole, or in part.
First published November 2011.

Scion and the Scion logo are trademarks or registered trademarks of Scion Corporation.

© 2011 Scion

Table of Contents

1.0 Revision History.....	5
2.0 Qualification Representative and Reviewer Details.....	6
2.1 Qualification Representative Details.....	6
2.2 Reviewer Details.....	7
2.3 Quality Assurance/Control Details.....	8
3.0 Customer Responsibilities.....	9
4.0 Qualification Guidelines and GMP Documentation.....	10
4.1 Qualification Summary.....	10
4.2 Qualification Guidelines.....	10
4.3 Page Numbering of Appendices.....	11
4.4 Exception Reports.....	12
4.5 Reference Documents.....	13
4.6 Required Materials.....	14
4.7 General Guidelines.....	15
4.8 Specific Instructions for Documentation.....	15
4.9 Documentation Corrections.....	16
4.10 Marking Procedures Not Applicable.....	16
4.11 Addendums.....	17
4.12 Addendum Example.....	17
5.0 Operational Qualification.....	18
5.1 OQ Preparation.....	18
5.2 System Description.....	20
5.3 Data Sheet Specifications.....	21
5.4 EI Precision Test TQ.....	22
5.5 EI Precision Test SQ.....	23
5.6 Final Evaluation.....	25
6.0 Protocol Approval.....	26
6.1 Protocol Acceptance / Approval by Customer.....	26

6.2 Operational Qualification Protocol Assignment 6.3	27
Protocol Acceptance / Protocol Approval by Scion 6.4	27
Remarks	28
Appendices	29
A.1 Qualification Representative Details	30
B.1 Exceptions	31
Operational Qualification Protocol Certification	36

1.0 Revision History

This qualification protocol is updated as necessary, which includes the event of any regulatory changes to Title 21 of the Code of Federal Regulations (21 CFR) Parts 210 and 211 (if applicable), any software or hardware changes, or updates that may impact on regulatory compliance.

Issue Number	Date	Comments

2.0 Qualification Representative and Reviewer Details


2.1 Qualification Representative Details

Each person responsible for executing any part of this Protocol must complete the table below, providing a sample of their signature and initials, and recording the date the Qualification was performed.

Qualification representatives are nominated to execute and verify the completeness of the test protocol and correctness of all entries.

All testing must be performed in accordance with procedures outlined in this manual. The representative must be trained and qualified to perform the procedures outlined in this document.

A copy of their appropriate qualifications is to be inserted into "Qualification Representative Details" on page 30.

Name (Print)	SOMCHAI POHTONGKAM
Title	ENGINEER
Signature	
Initials	SOMCHAI
Date	16 July 24

Name (Print)	
Title	
Signature	
Initials	
Date	

2.2 Reviewer Details

Each representative responsible for reviewing any part of this protocol must record their details in the following tables, providing a sample of their signature and initials, and recording the date the qualification was performed.

An employee or designee of the company operating the instrument must review these qualification procedures. All calculations and data will be checked by the reviewer. Data review must be performed in accordance with the qualification guidelines "Qualification Guidelines and GMP Documentation" on page 10 and in compliance with current Good Manufacturing Practice (cGMP) as specified by 21 CFR Parts 210 and 211.

Documentation supporting training in the area of data review and cGMP must be carefully maintained and reviewed by the instrument owner.

Reviewer representatives are responsible for reviewing the completeness of the qualification protocol and accuracy of all entries.

Name (Print)	
Title	
Signature	
Initials	
Date	

Name (Print)	
Title	
Signature	
Initials	
Date	

2.3 Quality Assurance/Control Details

As Quality Assurance/Control (QA/QC), who is empowered to approve instrument compliance documents, I approve the procedures in the SCION Operational Qualification Protocol, which I may have amended, I accept the qualification of the Qualification Representative, and I will review and initial the results.

Name (Print)	
Title	
Signature	
Initials	
Date	

Name (Print)	
Title	
Signature	
Initials	
Date	

3.0 Customer Responsibilities

The customer shall ensure that the Preventive Maintenance (PM) or Installation Qualification (IQ) up to point 9.11 is completed. A customer representative should be available at all times during the Operational Qualification Protocol.

Note The Operational Qualification Protocol test procedure should be performed after significant repairs, and at least once a year.

Qualification Rep. Initials	<i>Sachit P.</i>	Reviewer Initials		QA/QC Initials	
Date	<i>16 July 24</i>	Date		Date	

4.0 Qualification Guidelines and GMP Documentation

4.1 Qualification Summary

At the end of qualification execution, all tables and data entry fields must be completed and all test results, where specified, must be printed and attached to the protocol.

The Qualification Representative and the Reviewer must sign (signature or initials) and date each page that has a signature field. This represents agreement and acceptance of all data and information on the signed page.

Note Scion does not provide instructions for full Qualification of the personal computer (PC) used to operate the SCION. If further qualification of the PC is required the end-user must contact the PC manufacturer.

Note Scion does not provide full qualification instructions for non-Scion manufactured accessories. Limited instructions may be supplied. If qualification of a non-Scion accessory is required, the end user must contact the accessory manufacturer.

4.2 Qualification Guidelines

The following are general guidelines for performing the qualification tests in accordance with cGMP for the Manufacturing, Processing, Packaging, or Holding of Drugs per 21CFR Parts 210 and 211. Additional local requirements may also apply.

- Read the guidelines before starting the qualification.
- Perform all tests exactly as written.
- Use a pen with permanent blue or black ink unless otherwise specified by company policy.
- Neatly strike out any incorrect words or numbers, made while writing comments or recording results, information or data within this Protocol, with a single line. The word(s) crossed out must remain legible. Write the correction as close as possible to the original entry. Write a brief description of the error. For example, write 'Transcription error' or 'Re-written for clarity'. Initial and date the change.
- Entering initials where a signature is requested, and vice versa is permitted. The exception to this is in 2.0 : Qualification Representative and Reviewer Details on page 6, where examples of each person's signature and initials are required.
- Use the date format dd Mon yyyy (e.g. 08 Mar 2011) unless otherwise specified by company policy.

- Complete all tables and data fields to comply with this protocol. Blank fields are not permitted. For items that are not applicable, draw a line through the field, and write 'N/A' (Not Applicable). If entire tables or sections of tables are not applicable, strike a line either through the entire table or the specific area and enter 'N/A'. Complete the signature fields on the page.
- Write 'Pass', 'Fail' or 'N/A' as applicable to the test requirement or outcome.
- Ensure that results and/or specific documents are printed and attached to the specified appendix.
- The Qualification Representative and Reviewer must both sign (signature or initials) and date the signature fields on each page. This represents agreement and acceptance of all data and information on the page.

4.3 Page Numbering of Appendices

Each page that is inserted after the appendix is numbered with the letter of the appendix and a sequential number. The appendix page number must be initialed and dated by both the Qualification Representative and the Reviewer.

For example, pages inserted after Appendix C are numbered

C-1, C-2, C-3... etc. along with the initials and date.

If the reverse of each appendix page is left blank, it should be marked 'N/A' and signed and dated.

When the IQ is complete the total number of pages inserted after each appendix is written on the front page of the respective appendix sheet.

Qualification Rep. Initials	<i>Seah P.</i>	Reviewer Initials		QA/QC Initials	
Date	<i>16 July 24</i>	Date		Date	

4.4 Exception Reports

An exception to the protocol occurs when the observed result differs from the acceptance criteria or expected result.

All exceptions to the protocol must be documented in the Exception Report. The Exception Report includes a detailed description of the exception and resolution by the Qualification Representative.

Each Exception Report shall be issued with a unique identification number in the form ERID-XX-X. This number is generated by the page number on which the exception occurred followed by a sequential number indicating each exception found on the page.

For example, if an exception occurs on page 34, the Exception Report shall be identified as 'ERID-34-1'. If another exception occurs on page 34, the second report shall be identified as 'ERID-34-2'. This identification number should be recorded in the 'Pass / Fail / N/A' field after each test.

Each Exception Report must be signed by the Qualification Representative and the Reviewer as evidence of approval.

The Exception Report is inserted in the appropriately named appendix and numbered as per Section 4.3 of this protocol.

Qualification Rep. Initials	<i>Sadul P.</i>	Reviewer Initials		QA/QC Initials	
Date	<i>16 July 24</i>	Date		Date	

4.5 Reference Documents

The following documents are relevant to this Qualification:

- Installation Qualification Protocol
- Completed service report from Preventative Maintenance (PM) schedule

Qualification Rep. Initials	<i>S. P.</i>	Reviewer Initials		QA/QC Initials	
Date	<i>16 July 24</i>	Date		Date	

4.6 Required Materials

The following stock solutions are required:

- 100 fg/ μ L OFN 394204200
- 1 pg/ μ L OFN 393065201
- 100 pg/ μ L OFN 393110101
- 10 pg/ μ L BZP 93065301
- 100 pg/ μ L BZP 394206200

The above solutions will be used to prepare the following working solutions which will be required to execute this OQ:

Note Refer to Appendix 1 for the preparation of the standard solutions.

Qualification Rep. Initials	<i>Sachin P.</i>	Reviewer Initials		QA/QC Initials	
Date	<i>16 July 24</i>	Date		Date	

4.7 General Guidelines

The following are general cGMP guidelines.

- Perform each procedure exactly as written.
- Fill in each item on the form or mark it Not Applicable (N/A).
- If an item is marked N/A, initial it and date it.
- The Reviewer reviews and initials all entries recorded by the Qualification Representative.
- Keep all raw data. The Qualification Representative and the Reviewer will initial it, and date it.
- Do not destroy raw data.
- Attach raw data from an instrument, such as the SCION, as an Addendum using the instructions in the following Addendums section.
- If several instruments are qualified simultaneously, reference shared information, such as standard preparation and chemical information, to the document containing the original information by its model and instrument identification number.
- Label all reference standards as required by local regulations.
- Record the time each reference standard was opened.
- Use reference standards within 24 hours of preparation.

Qualification Rep. Initials	<i>Sabur P.</i>	Reviewer Initials		QA/QC Initials	
Date	16 July 24	Date		Date	

4.8 Specific Instructions for Documentation

The Reviewer designates specific documentation instructions as follows.

Permanent Ink Color	Blue
Preferred Date Format	16 July 24

If more instructions are required: Use an addendum sheet, write the addendum number, and a brief description.

Qualification Rep. Initials	<i>Sabur P.</i>	Reviewer Initials		QA/QC Initials	
Date	16 July 24	Date		Date	

4.9 Documentation Corrections

Note All original entries must remain legible after corrections are made.

1. Draw a line through the incorrect information.
2. Write the correction as close as possible to the original entry, or enter a footnote.
3. Write a brief description of the error. For example, write "transcription error," "rewritten for clarity," or "correcting wrong entry".
4. Initial and date the change.

Qualification Rep. Initials	<i>Sarah P.</i>	Reviewer Initials		QA/QC Initials	
Date	16 July 24	Date		Date	

4.10 Marking Procedures Not Applicable

Some sections may not be relevant for the qualification. To indicate that a procedure or part of a form is unnecessary and that it was not forgotten or inadvertently overlooked:

1. Draw a line through the portion that is not applicable. Write the letters N/A (for not applicable), your initials, and the date near the diagonal line.
2. If a procedural step is unnecessary, select N/A if it is indicated, or write a comment in an Addendum.

The Qualification Representative and the Reviewer enter their initials and the date near the line.

Note The Qualification Representative and Reviewer must sign and date all forms, even when part or all of the form is marked N/A.

Qualification Rep. Initials	<i>Sarah P.</i>	Reviewer Initials		QA/QC Initials	
Date	16 July 24	Date		Date	

4.11 Addendums

The following are reasons to complete an addendum sheet:

- A deviation needs documentation.
- Additional information or data needs to be recorded.
- Insufficient space to include the correction on the sheet where the error was made.

Qualification Rep. Initials	<i>Sarah P.</i>	Reviewer Initials		QA/QC Initials	
Date	<i>16 July 24</i>	Date		Date	

4.12 Addendum Example

The following is an example of using an addendum sheet to document a deviation.

If some of the items on the sales order were not present, you could do the following:

1. Use an addendum sheet.
2. Write Instrument Delivery on the Procedure line.
3. Write the addendum page number followed by a letter. For example: page 12A, where 12 is the page and A represents the first addendum on that page.
4. Write the plan to obtain the missing items, which may be the following:
 - Scion notified that Part Number XXXXX and XXXX are missing.
 - Scion replied that the parts are in stock and will be sent overnight. While
 - waiting for the parts to arrive, I will continue to set up the instrument.
5. Review the plan with the Reviewer and make the necessary modifications.
6. Document the arrival of the parts and write that this addendum is resolved. Attach a copy of delivery documents and create addendum pages as required.

Qualification Rep. Initials	<i>Sarah P.</i>	Reviewer Initials		QA/QC Initials	
Date	<i>16 July 24</i>	Date		Date	

5.0 Operational Qualification

This chapter contains the tests to be completed to perform an Operational Qualification for the SCION.

5.1 OQ Preparation

The following must be done before starting the OQ:

1. Preventative Maintenance must have been completed and signed off by the Qualification Representative, Reviewer, and QA/QC person, and attach a copy of the service report and add an addendum number.

Addendum P.M. Protocol

Qualification Rep. Initials	<i>Sachin P.</i>	Reviewer Initials		QA/QC Initials	
Date	16 July 24	Date		Date	

2. OQ must have been completed and signed off by the Qualification Representative, Reviewer, and QA/QC person.

Qualification Rep. Initials	<i>Sachin P.</i>	Reviewer Initials		QA/QC Initials	
Date	16 July 24	Date		Date	

3. The QA/QC person must review, approve, append (if necessary), and sign the Pre-execution Approval.

Qualification Rep. Initials	<i>Sachin P.</i>	Reviewer Initials		QA/QC Initials	
Date	16 July 24	Date		Date	

4. The Qualification Representative and the Reviewer must sign and date the Pre-execution Approval.

Qualification Rep. Initials	<i>S. P.</i>	Reviewer Initials		QA/QC Initials	
Date	<i>16 July 20</i>	Date		Date	

5.2 System Description

5.2.1 SCION Description

Installation Date: 2015	Principal Operator:	Phone Number:
Company Information		
Company Name: United Analyst and Engineering		Installation Site: LAB
Address: 3 Soi Udomsak 41		Address/Location: Sukhumvit Rd.
City, State: Bangkok, Prakhonong		City, State: Bangkok
Zip / Country: Thailand.		Zip / Country: 10260
System Description		
SCION SQ	Serial Number: GQS1203F021	
Sales Order Number:	Sales Order Addendum Number:	
GC		
Module Type: Scion A51	Serial Number: BR1203M099	
AutoSampler		
Module Type: CP8400	Serial Number: BR1203M331	
MS Workstation		
Version: MSWS 8-2.1	Serial Number: 01106-6711-BBQ-450d	
Computer Operating System		
Operating System: Windows 7	Version: Pro	Serial No.: 00386-150-436-158 Service Pack: -
Computer		
Make: Dell	Model: optiplex	Serial No.: DNNYHSI Hard Drive 1TB Size / RAM: 4GB
Addendum Number(s):	2. System description	

Qualification Rep. Initials	Santana P.	Reviewer Initials		QA/QC Initials	
Date	16 July 24	Date		Date	

5.3 Data Sheet Specifications

Run these tests after the instrument has pumped down and is leak free. Use the factory methods. Follow the Installation Procedure; complete this section and the appropriate line of the OQ Summary. Print out the methods and results and attach as addendums. Use the factory test column Br-5ms 15m X 250im X 0.25im.

Table 5-1 TQ Specification

Mode	Concentration	Scan Range	Result †	N/A	Pass	Fail	Addendum
EI Full Scan	1 pg OFN	50-300	S/N \geq 500:1				
EI MRM	100 fg OFN	272-222	S/N \geq 5000:1				
PCI Full Scan‡	10 pg BZP	80-230	S/N \geq 50:1			N/A	
NCI Full Scan‡	1 pg OFN	200-300	S/N \geq 4000:1				

† The Signal-to-Noise ratio S/N values are based on RMS noise figure.

‡ CI tests use methane gas as reagent gas.

For any tests that did not pass, complete an Addendum for each, write the Addendum number and a brief description.

Qualification Rep. Initials	<i>Sander P.</i>	Reviewer Initials		QA/QC Initials	
Date	<i>16 July 24</i>	Date		Date	

Table 5-2 SQ Specification

Mode	Concentration	Scan Range	Result †	N/A	Pass	Fail	Addendum
EI Full Scan	1 pg OFN	50-300	S/N \geq 600:1		✓		
PCI Full Scan‡	100 pg BZP	80-230	S/N \geq 600:1	✓			
NCI Full Scan‡	200 fg OFN	200-300	S/N \geq 1000:1	✓			

Qualification Rep. Initials	<i>Sander P.</i>	Reviewer Initials		QA/QC Initials	
Date	<i>16 July 24</i>	Date		Date	

5.4 EI Precision Test TQ

The following precision tests are for systems with autosamplers only. The test solution is 1 pg/ μ L OFN test mix part number 393065201.

The following is the required precision for 10 consecutive injections:

Injection	Retention Time	Peak Area
1		
2		
3		
4		
5		
6	N/A	
7		
8		
9		
10		
% RSD		

As an alternative, a % RSD summary report from MSWS can be added as an addendum.

Addendum N/A

	N/A	Pass	Fail	Addendum
Observed Mass is between 271.6 m/z to 272.4 m/z, which is ± 0.4 of the expected m/z.	✓			
Retention Time $\leq 1\%$ Relative Standard Deviation (RSD).	✓			
Peak Area $\leq 10\%$ Relative Standard Deviation (RSD).	✓			

To complete this section use the factory MRM method on the system CD. Print a copy of the method and add as an addendum.

Addendum N/A

If the hardware is not the same as the factory method, then note this in the addendum and how the hardware available has been configured to compensate. The most common variation here is the sampler, where the Combi Pal has been used instead of the 8400. This will have no impact on results and can be tracked and recorded in the addendum.

5.5 EI Precision Test SQ

The following precision tests are for systems with autosamplers only. The test solution is 1 pg/ μ L OFN test mix part number 393065201.

The following is the required precision for 10 consecutive injections:

Injection	Retention Time	Peak Area
1	3.645	39800
2	3.643	35585
3	3.640	38804
4	3.642	36679
5	3.641	35202
6	3.642	36225
7	3.641	37711
8	3.642	36452
9	3.642	35329
10	3.643	38639
% RSD	0.038	4.343

As an alternative, a % RSD summary report from MSWS can be added as an addendum.

Addendum _____

	N/A	Pass	Fail	Addendum
Observed Mass is between 271.6 m/z to 272.4 m/z, which is ± 0.4 of the expected m/z.		✓		
Retention Time $\leq 1\%$ Relative Standard Deviation (RSD).		✓		
Peak Area $\leq 10\%$ Relative Standard Deviation (RSD).		✓		

To complete this section use the factory Scan method on the system CD. Print a copy of the method and add as an addendum.

Addendum N/A

If the hardware is not the same as the factory method, then note this in the addendum and how the hardware available has been configured to compensate. The most common variation here is the sampler, where the Combi Pal has been used instead of the 8400. This will have no impact on results and can be tracked and recorded in the addendum.

Addendum N/A

5.6 Final Evaluation

	N/A	Pass	Fail	Addendum
Is the equipment in normal operation condition?		✓		
Have all of the OQ requirements been completed?		✓		

Qualification Rep. Initials	<i>Sudhal P.</i>	Reviewer Initials		QA/QC Initials	
Date	<i>16 July 24</i>	Date		Date	

6.0 Protocol Approval

6.1 Protocol Acceptance / Approval by Customer

I agree that the procedures and information referenced in this document are applicable.

Instrument(s):

Scion 451 SQ with CP8400

Serial Number(s):

GQS 1203F021

Sales Order Number:

Company Name:

United Analyst and Engineering Consultant Co., Ltd.

I agree that the Operational Qualification Protocol has been satisfactorily completed.



I confirm that the Operational Qualification Protocol has not been completed, because of these failed (non-passed) items



Authorized Customer Representative

Name (Print)	
Title	
Signature	
Initials	
Date	

6.2 Operational Qualification Protocol Assignment

This Operational Qualification Protocol document is used for:

Operational Qualification Protocol as final test at Scion	<input type="checkbox"/>
Operational Qualification Protocol after Installation Qualification	<input type="checkbox"/>
Operational Qualification Protocol after Preventive Maintenance and OQ completion.	<input checked="" type="checkbox"/>

6.3 Protocol Acceptance / Protocol Approval by Scion

I agree that the procedures and information referenced in this document are applicable.

Instrument(s): Scion 451 SQ with DP8400

Serial Number(s): GLQS1203F21

Sales Order Number:

Company Name: United Analyst and Engineering Consultant Co., Ltd.

Scion Certified Engineer

Name (Print)	SOMCHAI POMTONGKAM
Title	ENGINEER
Signature	<i>Somchai P.</i>
Initials	SOMCHAI
Date	

6.4 Remarks

Appendices

Each page that is inserted after the appendix is numbered with the letter of the appendix and a sequential number. The appendix page number must be initialed and dated by both the Qualification Representative and the Reviewer.

For example, pages inserted after Appendix C are numbered C-1, C-2, C-3... etc along with the initials and date.

If the reverse of each appendix page is left blank it should be marked NA and signed and dated.

When the OQ is complete the total number of pages inserted after each appendix is written on the front page of the respective appendix sheet.



Certificate

It is hereby certified that

Mr. Somchai Pohtongkam

Has successfully completed the Service & Application Training for

Scion Chromatography Products

Training Contents were:

**Hardware Operation, Software operation, Data analysis and Installation, &
Troubleshooting of Model:**

SCION GC, GCMS SQ, GCMS TQ

At Techcomp Singapore
By Mr. Michael Mei (Service Manager)

On 11th~15th July 2016

Hans van den Heuvel
Commercial Director
Scion Instruments

Date: 19 July 2016

Cert. No. : TSG-SCIONGC-15011602

เอกสารไม่ควบคุม

A.1 Qualification Representative Details

The Qualification Representative is to insert a copy of their appropriate qualification(s) after this page.

No. of Pages Inserted	
-----------------------	--

This area is intentionally left blank.

B.1 Exceptions

Each Exception Report shall be issued with a unique identification number in the form of ERID-XX-X. This number is generated by the page number on which the exception occurred followed by a sequential number indicating each exception found on the page.

For example, if an exception occurs on page 34, it shall be identified as Exception Report 'ERID-34-1'. If another exception occurs on page 34, the second exception shall be identified as 'ERID-34-2'. This identification number should be recorded in the pass/fail field after each test.

Insert Exception Reports (if any) after this page.

No. of Pages Inserted	N/A
-----------------------	-----

This area is intentionally left blank.

Addendum Procedure: P.M. Protocol Page Number: 1

Qualification Rep. Initials	<i>Salu' R</i>	Reviewer Initials		QA/QC Initials	
Date	<i>16 July 24</i>	Date		Date	

Addendum Procedure: 2. System description Page Number: 5

Qualification Rep. Initials	<i>Sauhal P.</i>	Reviewer Initials		QA/QC Initials	
Date	<i>16 July 24</i>	Date		Date	

Addendum Procedure: 3-Test Result Page Number: 30

Qualification Rep. Initials	<i>Sawden P.</i>	Reviewer Initials		QA/QC Initials	
Date	<i>16 July 24</i>	Date		Date	

Addendum Procedure: 4. Certificate Page Number: _____

Qualification Rep. Initials	<u>Sandhu P.</u>	Reviewer Initials		QA/QC Initials	
Date	<u>16 July 24</u>	Date		Date	

Operational Qualification Protocol Certification

for

SCION

with the serial number

GRS 1203 F 21

has successfully completed all criteria for hardware Operational Qualification Protocol
as detailed in this document.

Scion Certified Engineer

SOMCHAI POHTONGKAM

Name (please print)

Somchai P.

Signature

16 July 24

Date

Authorized Customer Representative

Name / Function (please print)

Signature

Date

Customer Address

United Analyst and Engineering Consultant co.,Ltd.