

ภาคผนวก ข

เอกสารสอบเทียบเครื่องมือ

List of Instrument Certificates for Environmental Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*
1	Analytical Balance	FAT OIL AND GREASE	Mettler Toledo	AB204-S/FACT / 1129361010	United Analyst and Engineering Consultant Co., Ltd.	250422 1 BL002 25	23/4/2025	22/4/2026
2	Analytical Balance	TOTAL DISSOLVED SOLIDS	Mettler Toledo	XSR205DU / C210685394	National Food Institute,Ministry of Industry, Thailand	2502226-002-01	20/3/2025	19/3/2026
3	Analytical Balance	TOTAL SUSPENDED SOLIDS	Mettler Toledo	XSR205DU / C009071872	National Food Institute,Ministry of Industry, Thailand	2502226-001-01	20/3/2025	19/3/2026
4	BOD Incubator	BIOCHEMICAL OXYGEN DEMAND	ARCO	UC4-1320 / -	Technology Promotion Association (Thailand-Japan)	25TM1003	8/7/2025	7/7/2026
5	DO Meter	BIOCHEMICAL OXYGEN DEMAND	YSI	5100 / 11B 101863	Technology Promotion Association (Thailand-Japan)	25TW29	17/2/2025	16/2/2026
6	Digestion Units	TOTAL KJELDAHL NITROGEN	Foss Tecator	2520 Auto / 91794469	National Food Institute Ministry of Industry, Thailand	2501440-001-01	27/1/2025	26/1/2026
7	Hot Air Oven	TOTAL DISSOLVED SOLIDS TOTAL SUSPENDED SOLIDS	Memmert	UF55 / B212.0411	Technology Promotion Association (Thailand-Japan)	25TM579	19/3/2025	18/3/2026
8	Kjeltec System Distilling Unit	TOTAL KJELDAHL NITROGEN	Foss Tecator (Labtec)	KT200 / 91790524	FOSS South East Asia	13319	27/1/2025	26/1/2026
9	Kjeltec Distillation Unit	TOTAL KJELDAHL NITROGEN	FOSS	Kjeltec 8100 / 91889052	FOSS South East Asia	13854	24/2/2025	23/2/2026
10	pH Meter	pH	Horiba	LAQUA-PH210 / HA1M0043	technology promotion association (thailand-japan	25CH263	28/2/2025	27/2/2026

Due Date of Calibration* : Based on the annual calibration plan. At least 1 time per year.

Certificate of Calibration

Certificate No.: 250422-1-BL002-25

Code No.: BL002-25

Page: 1 of 3

Customer Name: United Analyst and Engineering Consultant Co., Ltd.
Address: 3 Soi Udom suk 41, Sukhumvit Rd., Bang Chak, Phar Khanong, Bangkok 10260

Equipment: Electronic Balance

Manufacturer: Mettler Toledo

Model: AB204-S/FACT

Serial No.: 1129361010

Asset No.: UAE.WAS.002/2552

Building: N/A **Floor:** 1 **Room:** 107

Received Date: April 22, 2025

Date of Calibration: April 23, 2025

Calibration Conditions: Temperature 22.8 °C to 23.4 °C
Humidity 54.8 % to 68.9 %
Pressure 756.6 mmHg to 758.2 mmHg

Calibrated by: Sakkarin Srirahang

Approved by: Suwit Chotnok

Signature:

Issued Date: April 25, 2025

Note: 1) The Uncertainties are for a confidence probability of approximately 95%

2) This Certificate is valid only to the item calibrated on date and place of calibration.

3) This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the United Analyst and Engineering Consultant Co., Ltd. (UAE)

เอกสารไม่ควบคุม

Certificate No.: 250422-1-BL002-25

Code No.: BL002-25

Page: 3 of 3

Equipment: Electronic Balance **Manufacturer:** Mettler Toledo
Model: AB204-S/FACT **Readability:** 0.0001 g
Serial No.: 1129361010 **ID No.:** UAE.WAS.002/2552
Max. Capacity: 220 g
Calibration Date: April 23, 2025

Calibration Result: (Continued)

Calibration Range: 0 - 200 g

Calibration Adjustment: Internal Calibration

3. Error of indication from nominal or conventional mass value:

Nominal Value (g)	Reference Value (g)	Indication (g)	Correction (g)	Uncertainty (\pm mg)	Coverage Factor k
Unload	0.0000000	0.0000	0.0000	0.10	2.05
0.01	0.0100025	0.0099	0.0001	0.10	2.05
0.05	0.0500056	0.0500	0.0000	0.10	2.05
0.1	0.1000012	0.0999	0.0001	0.10	2.05
0.5	0.5000133	0.5000	0.0000	0.10	2.05
1	1.0000105	1.0000	0.0000	0.10	2.05
10	10.000010	10.0000	0.0000	0.11	2.04
40	40.000076	40.0000	0.0000	0.14	2.00
50	50.000056	50.0000	0.0001	0.13	2.00
80	80.000107	80.0000	0.0001	0.18	2.00
100	100.000109	99.9999	0.0002	0.17	2.00
120	120.00015	119.9999	0.0003	0.21	2.00
150	150.000165	149.9998	0.0003	0.24	2.00
160	160.000175	159.9997	0.0005	0.26	2.00
200	200.000129	199.9998	0.0004	0.30	2.00

4. Effect of Tare test:

Tare Load (g)	Test Load (g)	Indication (g)	Correction (g)
100	20.000041	19.9999	0.0001
	40.000076	39.9998	0.0002
	60.000066	59.9997	0.0003
	80.000107	79.9999	0.0002
	100.000168	100.0004	-0.0003

Remark:

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k, providing a confidence probability of approximately 95%.

o---o-End-o---

เอกสารไม่ควบคุม

Certificate No.: 250422-1-BL002-25

Code No.: BL002-25

Page: 2 of 3

Equipment: Electronic Balance **Manufacturer:** Mettler Toledo
Model: AB204-S/FACT **Readability:** 0.0001 g
Serial No.: 1129361010 **ID No.:** UAE.WAS.002/2552
Max. Capacity: 220 g
Calibration Date: April 23, 2025
Condition As-Received: In Condition

Condition of Equipment:

Condition of This Result of Calibration:

1. Calibration Method: This instrument was calibrated by method UAE.CP.CAL.006 In-House Method based on UKAS Lab 14 : 2022

2. Reference Standards:

Reference Standard:	Model	Serial No.	Calibrated By	Certificate No.	Traceability	Due Date
Standard Weight Class E2 (OIML)	1 mg to 1 kg	8749109122	AMARC	25-009359	Mettler-Toledo	21-Jan-27
Standard Weight Class F1 (OIML)	1 mg to 200 g	11119512	AMARC	24-013840	Mettler-Toledo	04-Feb-26
Instrument	Model	Serial No.	Calibrated By	Certificate No.	Traceability	Due Date
Thermo-Hygro-Baro Meter	MHB-3825D	AK.46457	SUCCESS	SG-H-00997/67	Success Gateway	21-Nov-25
Thermo-Hygro-Baro Meter	MHB-3825D	AK.46457	TPA	25P795	TPA	25-Feb-26

3. This certification is traceable to SI Unit

4. This certification was certified only for the instrument we calibrated

5. This result of calibration was found accurate as shown on date and place of calibration only.

6. Through the reference standard laboratory of AMARC 25-009359 Calibration 0152

Calibration Result:


1. Repeatability of Reading:

Nominal Value (g)	Standard Deviation of Reading (g)
200*	0.000045

2. Eccentric or off-center loading

A mass of 100 g was placed and moved to various position on pan

The Balance reading obtained is given in the table.

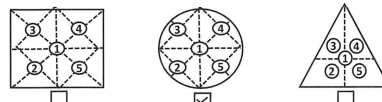

1 (g)	2 (g)	3 (g)	4 (g)	5 (g)	Maximum Difference (g)
100.0000	99.9996	99.9997	100.0003	100.0005	0.0005

เอกสารไม่ควบคุม

Calibration Certificate

Certificate No.: 2502226-002-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260

Page 1 of 4

Equipment: Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Serial No.: C2106085394

ID No.: UAE.WAO.010/2565

Order No.: 2502226

Operation No.: 2502226-002

Date of Receipt: 19 March 2025

Date of Calibration: 20 March 2025

Calibrated by Mr.Yothin Charoensuk
Scientist

Approved by (Mr.Pheraphat Tuanjit)
Manager, Division of Calibration Laboratory
Responsible for the Technical Management Team

Date of Issue: 25 March 2025

The uncertainties are for a confidence probability of approximately 95%

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

E-CO-009 Revision: 01 Date: 20-04-65

2008 ต.บางนา-สุขุมวิท 36 แขวงบางนา-สุขุมวิท เขตบางนา กรุงเทพมหานคร 10260
2008 Soi 36, Arun Amarin Road, Bang Yi Khan Subdistrict, Bang Phat District, Bangkok 10260, Thailand
Tel :+66(0) 2422 6688 Fax :+66(0) 2422 6545



nfi.or.th

เอกสารไม่ควบคุม

Calibration Report

Certificate No.: 2502226-002-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C210685394

ID No.: UAE.WAO.010/2565

Capacity: 82 g / 220 g

Date of Calibration: 20 March 2025

Page 2 of 4

Environment Condition: Ambient Temperature: 21.2 ± 0.6 °C Relative Humidity: 48 ± 3.5 %

Place of Calibration: 208 Balance Room, UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.

Condition of Equipment: Good Condition

Condition of This Results of Calibration:

1. Calibration Method: NFI Method W-MA-001 In-House Method based on UKAS Lab 14 : 2019

2. Reference Standards:

Reference Standard	Model	Serial No.	Calibrated By	Certificate No.	Due Date
Standard Weight Class E2	1mg to 200g	B505567572	TCS	M24041005	19 April 2025
Instrument	Model	Serial No.	Calibrated By	Certificate No.	Due Date
Thermo-Hygro Meter	608-H1	NF1.BTH 017/23	Quality Reborn	QR25-0542	10 February 2026

3. This certification is traceable to SI UNIT

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

Calibration Results:


1. Repeatability of Reading:

Nominal Value (g)	Standard Deviation of Reading (g)
40	0.0000042
80	0.0000042
100	0.0000000
200	0.0000000

2. Off-Center Error:

A mass of 100 g was placed and moved to various position on pan.

The balance reading obtained is given in the table.

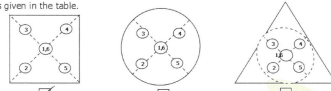

1	2	3	4	5	6	(Maximum Difference)
(g)	(g)	(g)	(g)	(g)	(g)	(g)
100.0001	100.0001	100.0001	100.0001	100.0001	100.0001	0.0000

F-CS-012 Revision: 01 Date: 20-04-65

2008 3668 ถนนสุขุมวิท 36 แขวงคลองเตย เขตวัฒนา กรุงเทพมหานคร 10110
2008 Soi 36, Arun Amarin Road, Bang Yi Khan Subdistrict, Bang Phlat District, Bangkok 10700, Thailand
Tel :+66(0) 2422 8668 Fax :+66(0) 2422 8545



nfi.or.th

Calibration Report

Certificate No.: 2502226-002-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C210685394

ID No.: UAE.WAO.010/2565

Capacity: 82 g / 220 g

Date of Calibration: 20 March 2025

Page 3 of 4

Calibration Results: (Continued)

Calibration Range: 0-80 g

Calibration Adjustment: Internal Calibration

3. Departure from Nominal Value: (Range: 0 - 82 g ; Resolution: 0.00001 g)

Nominal Value (g)	Standard Value (g)	Average Reading (g)	Correction (g)	Uncertainty (± g)	Coverage Factor k
Unload	0.000000	0.000000	0.000000	0.0000087	2.00
0.001	0.001003	0.001000	0.000000	0.0000090	2.00
0.005	0.005002	0.005001	-0.000001	0.0000092	2.00
0.01	0.010003	0.010002	-0.000002	0.0000089	2.00
0.05	0.049996	0.050001	-0.000001	0.0000096	2.00
0.1	0.100011	0.100002	-0.000001	0.0000111	2.00
0.5	0.500016	0.500004	-0.000002	0.000014	2.00
1	1.000003	1.000005	-0.000005	0.000016	2.00
2	2.000023	2.000006	-0.000004	0.000017	2.00
5	5.000015	5.000006	-0.000005	0.000020	2.00
10	10.000009	10.000005	-0.000004	0.000026	2.00
20	20.000030	20.000007	-0.000004	0.000037	2.00
30	30.000039	30.000009	-0.000005	0.000050	2.00
50	50.000028	50.000008	-0.000005	0.000068	2.00
80	80.000067	80.000013	-0.000006	0.000111	2.00

F-CS-012 Revision: 01 Date: 20-04-65

2008 3668 ถนนสุขุมวิท 36 แขวงคลองเตย เขตวัฒนา กรุงเทพมหานคร 10110
2008 Soi 36, Arun Amarin Road, Bang Yi Khan Subdistrict, Bang Phlat District, Bangkok 10700, Thailand
Tel :+66(0) 2422 8668 Fax :+66(0) 2422 8545



nfi.or.th

Calibration Report

Certificate No.: 2502226-002-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C210685394

ID No.: UAE.WAO.010/2565

Capacity: 82 g / 220 g

Date of Calibration: 20 March 2025

Page 4 of 4

Calibration Results: (Continued)

Calibration Range: >80-200 g

Calibration Adjustment: Internal Calibration

3. Departure from Nominal Value: (Range: >80 - 200 g ; Resolution: 0.0001 g)

Nominal Value (g)	Standard Value (g)	Average Reading (g)	Correction (g)	Uncertainty (± g)	Coverage Factor k
90	90.000010	90.00002	-0.000001	0.000015	2.00
100	100.000006	100.00001	0.000000	0.000016	2.00
110	110.000007	110.00002	-0.000001	0.000017	2.00
120	120.000009	120.00002	-0.000001	0.000018	2.00
130	130.000010	130.00002	-0.000001	0.000019	2.00
140	140.000013	140.00002	-0.000001	0.000019	2.00
150	150.000009	150.00002	-0.000001	0.000021	2.00
160	160.000010	160.00002	-0.000001	0.000022	2.00
170	170.000012	170.00002	-0.000001	0.000023	2.00
200	200.000013	200.00002	-0.000001	0.000028	2.00

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k , providing a level of confidence of approximately 95 %.

----- End -----

F-CS-012 Revision: 01 Date: 20-04-65

2008 3668 ถนนสุขุมวิท 36 แขวงคลองเตย เขตวัฒนา กรุงเทพมหานคร 10110
2008 Soi 36, Arun Amarin Road, Bang Yi Khan Subdistrict, Bang Phlat District, Bangkok 10700, Thailand
Tel :+66(0) 2422 8668 Fax :+66(0) 2422 8545



nfi.or.th

Calibration Certificate

Certificate No.: 2502226-001-01

Client name:

UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.

Address:

3 Soi Udomsuk 41, Sukhumvit Road,
Bangchack, Prakhonong, Bangkok 10260

Equipment:

Electronic Balance

Manufacturer:

METTLER TOLEDO

Model:

XSR205DU

Serial No.:

C009071872

ID No.:

UAE.WAO.012/2563

Order No.:

2502226

Operation No.:

2502226-001

Date of Receipt:

19 March 2025

Date of Calibration:

20 March 2025

Calibrated by Mr.Yothin Charoensuk
Scientist

Approved by (Mr.Pheraphat Tuanjit)

Manager, Division of Calibration Laboratory

Date of Issue: 25 March 2025

Responsible for the Technical Management Team

The uncertainties are for a confidence probability of approximately 95%

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65

2008 3668 ถนนสุขุมวิท 36 แขวงคลองเตย เขตวัฒนา กรุงเทพมหานคร 10110
2008 Soi 36, Arun Amarin Road, Bang Yi Khan Subdistrict, Bang Phlat District, Bangkok 10700, Thailand
Tel :+66(0) 2422 8668 Fax :+66(0) 2422 8545



nfi.or.th

Calibration Report

Certificate No.: 2502226-001-01

Equipment: Electronic Balance
Model: XSR205DU
Serial No.: C009071872
Capacity: 82 g / 220 g

Manufacturer: METTLER TOLEDO
Resolution: 0.00001 g / 0.0001 g
ID No.: UAE.WAO.012/2563

Date of Calibration: 20 March 2025

Page 2 of 4

Environment Condition: Ambient Temperature: 21.2 ± 0.6 °C Relative Humidity: 48 ± 3.5 %
Place of Calibration: 208 Balance Room, UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.

Condition of Equipment: Good Condition

Condition of This Results of Calibration:

1. Calibration Method: NFI Method W-MA-001 In-House Method based on UKAS Lab 14 : 2019

2. Reference Standards:

Reference Standard	Model	Serial No.	Calibrated By	Certificate No.	Due Date
Standard Weight Class E2	1mg to 200g	B505567572	TCS	M24041005	19 April 2025
Instrument	Model	Serial No.	Calibrated By	Certificate No.	Due Date
Thermo-Hygro Meter	608-H1	NFI.BTH 017/23	Quality Reborn	QR25-0542	10 February 2026

3. This certification is traceable to SI UNIT

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

Calibration Results:

1. Repeatability of Reading:

Nominal Value (g)	Standard Deviation of Reading (g)
40	0.0000052
80	0.0000042
100	0.0000000
200	0.0000000

2. Off-Center Error:

A mass of 100 g was placed and moved to various position on pan.

The balance reading obtained is given in the table.

1	2	3	4	5	6
(g)	(g)	(g)	(g)	(g)	(g)
100.0001	100.0001	100.0001	100.0001	100.0001	100.0002
					(Maximum Difference)
					0.0001

F-CS-012 Revision: 01 Date: 20-04-65

2008 ปตท. จำกัด (มหาชน) 36 หมู่ 5 ต.บางนา-สวนหลวง กรุงเทพมหานคร 10700
2008 Soi 36, Arun Amarin Road, Bang Yi Khan Subdistrict, Bang Phai District, Bangkok 10700, Thailand
Tel: +66(0) 2422 8688 Fax: +66(0) 2422 8545



nfi.eath

Calibration Report

Certificate No.: 2502226-001-01

Equipment: Electronic Balance
Model: XSR205DU
Serial No.: C009071872
Capacity: 82 g / 220 g

Manufacturer: METTLER TOLEDO
Resolution: 0.00001 g / 0.0001 g
ID No.: UAE.WAO.012/2563

Date of Calibration: 20 March 2025

Page 3 of 4

Calibration Results: (Continued)

Calibration Range: 0-80 g

Calibration Adjustment: Internal Calibration

3. Departure from Nominal Value: (Range: 0 - 82 g ; Resolution: 0.00001 g)

Nominal Value (g)	Standard Value (g)	Average Reading (g)	Correction (g)	Uncertainty (± g)	Coverage Factor k
unloaded	0.000000	0.000000	0.000000	0.0000089	2.00
0.001	0.001003	0.001000	0.000000	0.0000092	2.00
0.005	0.005002	0.005000	0.000000	0.0000094	2.00
0.01	0.010003	0.010000	0.000000	0.0000091	2.00
0.05	0.049996	0.050000	0.000000	0.0000098	2.00
0.1	0.100011	0.100000	0.000011	0.000011	2.00
0.5	0.500016	0.500000	0.000016	0.000014	2.00
1	1.000003	1.000001	-0.000001	0.000016	2.00
2	2.000023	2.000005	-0.000018	0.000017	2.00
5	5.000015	5.000005	-0.000010	0.000021	2.00
10	10.000009	10.000005	-0.000004	0.000026	2.00
20	20.000030	20.000012	-0.000018	0.000037	2.00
30	30.000039	30.000012	-0.000027	0.000050	2.00
50	50.000028	50.000014	-0.000014	0.000068	2.00
80	80.000067	80.000020	-0.000047	0.00011	2.00

F-CS-012 Revision: 01 Date: 20-04-65

2008 ปตท. จำกัด (มหาชน) 36 หมู่ 5 ต.บางนา-สวนหลวง กรุงเทพมหานคร 10700
2008 Soi 36, Arun Amarin Road, Bang Yi Khan Subdistrict, Bang Phai District, Bangkok 10700, Thailand
Tel: +66(0) 2422 8688 Fax: +66(0) 2422 8545



nfi.eath

Calibration Report

Certificate No.: 2502226-001-01

Equipment: Electronic Balance
Model: XSR205DU
Serial No.: C009071872
Capacity: 82 g / 220 g

Manufacturer: METTLER TOLEDO
Resolution: 0.00001 g / 0.0001 g
ID No.: UAE.WAO.012/2563

Date of Calibration: 20 March 2025

Page 4 of 4

Calibration Results: (Continued)

Calibration Range: >80-200 g

Calibration Adjustment: Internal Calibration

3. Departure from Nominal Value: (Range: >80 - 200 g ; Resolution: 0.0001 g)

Nominal Value (g)	Standard Value (g)	Average Reading (g)	Correction (g)	Uncertainty (± g)	Coverage Factor k
90	90.00010	90.00002	-0.00008	0.00015	2.00
100	100.00005	100.00001	-0.00004	0.00016	2.00
110	110.00007	110.00001	-0.00006	0.00017	2.00
120	120.00009	120.00002	-0.00007	0.00018	2.00
130	130.00010	130.00002	-0.00008	0.00019	2.00
140	140.00013	140.00002	-0.00011	0.00019	2.00
150	150.00009	150.00002	-0.00007	0.00021	2.00
160	160.00010	160.00002	-0.00008	0.00022	2.00
170	170.00012	170.00002	-0.00010	0.00023	2.00
200	200.00013	200.00002	-0.00011	0.00028	2.00

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k , providing a level of confidence of approximately 95 %.

----- End -----

F-CS-012 Revision: 01 Date: 20-04-65

2008 ปตท. จำกัด (มหาชน) 36 หมู่ 5 ต.บางนา-สวนหลวง กรุงเทพมหานคร 10700
2008 Soi 36, Arun Amarin Road, Bang Yi Khan Subdistrict, Bang Phai District, Bangkok 10700, Thailand
Tel: +66(0) 2422 8688 Fax: +66(0) 2422 8545



nfi.eath

Certificate of Calibration

Cert. No.: 25TM1003
Page : 1 of 3

Equipment : BOD Incubator

Manufacturer : ARCO

Model : UC4-1320

Serial No. : -

ID No. : UAE.WAO.026/2566

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260

Location : Lab Floor 1

Received Order : 07 July 2025

Calibration Date : 08 July 2025

Ambient Temperature : (26 ± 10) °C

Relative Humidity : (50 ± 30) %

AC Line Voltage : (220 ± 22) V

Calibrated by : Man Pattanapongpaiboon

Approved by :

Approved Signatory

() Chakrit Waewwanjua

() Suwit Imjai

(✓) Kunchit Promprat

Issue Date : 17 July 2025

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2507-0146OC-4
Procedure Used :-

Cert. No.: 25TM1003
Page : 2 of 3

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD).
The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Data Acquisition	MY59003411	24LM192	TPA	24 Dec 2025

2. This certificate is valid only to the item calibrated on date and place of calibration.


3. This measurement result is traceable to the International System of Unit maintained through :

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Not Available

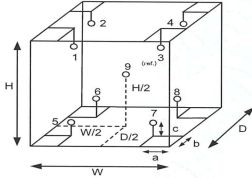

Environment during calibration		
	Beginning	Finished
Temp. (°C)	26	26
REL.Humid. (%)	55	57
AC Supply (Volt)	223	226

Position :	Ref. Std. ID No.:
1	25-20RTD-2/1
2	25-20RTD-2/2
3	25-20RTD-2/3
4	20RTD-2/4
5	20RTD-2/5
6	20RTD-2/6
7	20RTD-2/7
8	20RTD-2/8
9 (ref.)	20RTD-2/9

Probe Installation Details :

a = 10 cm
b = 10 cm
c = 10 cm

Dimension of Chamber :

D = 0.62 m
W = 1.2 m
H = 1.2 m
Capacity = 0.89 m³

เอกสารไม่ควบคุม



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2507-0146OC-4
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Not Available

Cert. No.: 25TM1003
Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor k
20.0	20.0	19.8	0.51	0.31	1.2	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	1	2	3	4	5	6	7	8	9 (ref.)	
20.0	19.958	19.896	20.018	20.116	19.830	19.771	19.871	19.880	19.909	0.75

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3 : EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL. 0-2717-3000 FAX. 0-2719-9484

Certificate of Testing

Cert.No.: 25TW29
Page.: 1 of 2

Equipment : DO Meter
Manufacturer : YSI
Model : 5100
Serial No. : 11B 101863
ID No. : UAE.WAO.004/2554
Received Date : 14 February 2025
Test Date : 17 February 2025
Reference : 2502-0473DSC-1
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260

Laboratory Condition : Temperature (25 ± 5) °C
Humidity (50 ± 20) %
Test Procedure : In - house method : CP-CH9
by Comparison Technique with Azide Modification Method

Tested by : Walalak Sirinthean

Approved by : [Signature]

Approved Signatory

() Chakrit Waewwanjua
() Ponpan Paipim
(✓) Saithip Meangmai

Issue Date : 18 February 2025

เอกสารไม่ควบคุม



Cert.No.: 25TW29
Page.: 2 of 2

Condition of this result of calibration

1. Reference Standard Instruments :

This certification is traceable to the International System of Unit through the reference standards laboratory of Industrial Calibration Center, Technology Promotion Association (Thailand-Japan).

Instruments	Serial No.	ID No.	Certificate No.	Due Date
1. Burette	-	130BU10	23CG1172	22 Mar 2025
2. Balance	14233821	110RC001	24MM131	04 July 2025

2. Standard Material :-

Material	Manufacturer	Lot.No.	Assay
Sodium Thiosulfate 5-Hydrate AR	KEMAUS	2203162447	99.6%

Result : Dissolved Oxygen Meter Adjustment With Air 100 %
Dissolved Oxygen Probe No.: 24F100202

Titration Method (Azide Modification Method) (mg/L)	DO Meter Reading (mg/L)	Standard Deviation (mg/L)
8.22	8.22	0.0055

This report was certified only for the instrument we tested. It is allowable to use for study intend to use for advertising and referral purpose is prohibited. This report may not be reproduced other in full, without written approval of the laboratory

-o0o-

เอกสารไม่ควบคุม

Verification Certificate

Certificate No.: 2501440-001-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road,
 Bangchack, Prakanong, Bangkok 10260

Page 1 of 4

Equipment:	Digestion Unit (Heating Block)
Manufacturer:	FOSS
Model:	Tecator Digestor 2520
Serial No.:	91905060
ID No.:	UAE.WAS.030/2566
Order No.:	2501440
Operation No.:	2501440-001
Date of Receipt:	27 January 2025
Date of Calibration:	27 January 2025

Calibrated by Mr. Worapob Sookthong
Scientist

Approved by (Mr. Phraphat Tuanjit)
Manager, Division of Calibration Laboratory
Responsible for the Technical Management Team

Date of Issue: 29 January 2025

The uncertainties are for a confidence probability of approximately 95 %.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65

เอกสารไม่ควบคุม

เอกสารไม่ครบ



Verification Report

Certificate No.:	2501440-001-01
Equipment:	Digestion Unit (Heating Block) Model: Tecator Digester 2520 Serial No.: 91905060 Resolution: 1 °C ID No : UAE.WAS.030/2566 Manufacturer: FOSS
Date of Calibration:	27 January 2025
Location:	Dry Laboratory (312), UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Environment Condition:	Ambient Temperature (29 ± 1) °C Relative Humidity (58 ± 2) % Line Voltage (224 ± 1) Volt

Condition of this results of Calibration:

1. This instrument was calibrated by insert standard thermocouples type R into its Digestion blocks and Calibration according to NFI Method W-TE-026 based on BS 4309 : 1968
- The Temperature scale used was based on ITS - 90 .
 - All data show below were final values and the initial data may be obtained upon request.

2. Reference Standard Instrument :

Instrument	Model	Serial No.	Certificate No.	Due Date	Through
Digital Thermometer with Thermocouple	34970A	MY44045576/MY41194453	TC24/0063	5-Jun-2025	N.M. Technical Center Laboratory
	Type R	S/CH1, R/CH2, R/CH3			

3. This certificate is traceable to International System of Units (SI Units).
4. This certificate was certified only for the instrument we calibrated.
5. This result of calibration was found accurate as shown on date and place of calibration only.
6. Condition of Calibrated item : Good

UUC* Description							
Time of Record	1	Hour	6	Minute	At	380	°C

7. Result of Calibration : ☒ Without adjustment ☐ After adjustment

F-CS-012 Revision: 01 Date: 20-04-65

UNCF M.6.4.D15-1(0)04-JUL-2023



Verification Report

Certificate No.: 2501440-001-01

Equipment: Digestion Unit (Heating Block)

Model: Tecator Digester 2520 Serial No.: 91905060

Resolution: 1 °C ID No : UAE.WAS.030/2566

Manufacturer: FOSS

Page 3 of 4

Date of Calibration:	27 January
Calibration point:	380 °C

Table 1 : Reporting of Temperature

Block No.	UUC* Setting (°C)	UUC* Reading (°C)	Stability (±°C)	Standard Thermometer (°C)	Uncertainty (±°C)
1	380	380	0.22	377.84	2.0
2	380	380	0.19	378.68	2.0
3	380	380	0.13	378.70	2.0
4	380	380	0.12	379.82	2.0
5	380	380	0.20	381.01	2.0
6	380	380	0.16	380.48	2.0
7	380	380	0.16	378.22	2.0
8	380	380	0.19	377.99	2.0
9	380	380	0.09	378.48	2.0
10	380	380	0.15	378.17	2.0
11	380	380	0.18	377.64	2.0
12	380	380	0.11	379.27	2.0
13	380	380	0.13	378.14	2.0
14	380	380	0.25	379.11	2.0
15	380	380	0.15	379.83	2.0
16	380	380	0.18	378.05	2.0
17	380	380	0.31	378.44	2.0
18	380	380	0.18	378.29	2.0
19	380	380	0.17	378.41	2.0
20	380	380	0.13	379.24	2.0

Note:

- UUC* = Unit Under Calibration
- Immersion depth of standard thermometer in tube level high of sand is equal heater plate of UUC.
- Stability = One-half of the greatest maximum difference of measured temperatures at one sensors, for at least half an hour after reaching steady state.

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor $k=2$ providing a level of confidence of approximately 95 %.

E-CS-012 Revision: 01 Date: 20-04-65


เอกสารไม่ครบคุณ



Verification Report

Certificate No.: 2501440-001-01
Equipment: Digestion Unit (Heating Block)
Model: Tecator Digestor ; Serial No.: 91905060
Resolution: 1 °C ID No.: UAE.WAS.030/2566
Manufacturer: FOSS
Date of Calibration: 27 January 2025 Page 4 of 4
Calibration point: 380 °C
Calibration result: Continued

Figure 1. Location of Reference Standard and Block Diagram of Digestion Unit

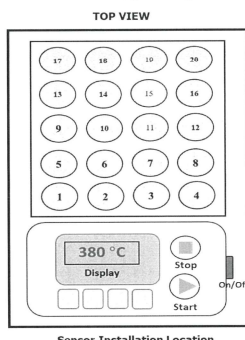

Sensor Installation Location

----- End -----

FCS-012 Revision: 01 Date: 20-04-65

เอกสารไม่ควบคุม

Certificate of Calibration

Cert. No.: 25TM579
Page : 1 of 3

Equipment : Hot Air Oven
Manufacturer : Memmert
Model : UF 55
Serial No. : B212.0411
ID No. : UAE.WAO.005/2556
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260
Location : Lab Floor 2
Received Order : 19 March 2025
Calibration Date : 19 March 2025
Ambient Temperature : (26 ± 10) °C
Relative Humidity : (50 ± 30) %
AC Line Voltage : (220 ± 22) V
Calibrated by : Man Pattanapongpaiboon
Approved by :
() Chakrit Waewwanjua
() Suwit Imjai
(✓) Kunchit Promprat
Issue Date : 27 March 2025

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : Hot Air Oven
Condition As-Received : Used Item
Reference : 2503-0437OC-3
Procedure Used :-
Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD) and Thermocouple Type T.

Cert. No.: 25TM579
Page : 2 of 3

The temperature scale used was based on ITS-90.

Condition of this result of calibration


- Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Data Acquisition	MY44073381	24LM73	TPA	18 May 2025
- This certificate is valid only to the item calibrated on date and place of calibration.
- This certification is traceable to the International System of Unit.

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Close

Environment during calibration		
	Beginning	Finished
Temp. (°C)	27	28
REL Humid. (%)	49	55
AC Supply (Volt)	221	224

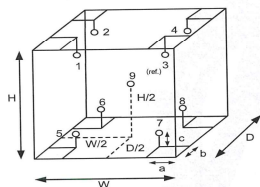

Probe Installation Details :
a = 5.0 cm
b = 5.0 cm
c = 5.0 cm
Dimension of Chamber :
D = 0.50 m
W = 0.80 m
H = 0.75 m
Capacity = 0.30 m³

Ref. Std. ID No.: @ Calibration Point		
Position :	(120,180) °C	(104) °C
1	23-01TC-01	1RTD-2/1
2	23-01TC-02	1RTD-2/2
3	23-01TC-03	22-01RTD-03
4	23-01TC-04	1RTD-2/4
5	23-01TC-05	1RTD-2/5
6	23-01TC-06	1RTD-2/6
7	23-01TC-07	23-01RTD-07
8	23-01TC-08	1RTD-2/8
9 (ref.)	23-01TC-09	23-01RTD-09



Equipment : Hot Air Oven
Condition As-Received : Used Item
Reference : 2503-0437OC-3
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Close

Cert. No.: 25TM579
Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor k
104.0	104.0	104.0	0.040	0.43	0.78	2
120.0	120.0	120.0	0.64	1.3	1.6	2
180.0	180.0	180.0	0.49	1.5	1.8	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	1	2	3	4	5	6	7	8	9 (ref.)	
104.0	104.335	104.135	104.363	104.317	103.649	103.738	104.179	104.229	104.025	0.42
120.0	119.575	119.366	119.807	119.905	118.994	119.194	119.888	119.994	120.064	1.1
180.0	180.286	179.510	180.401	180.551	179.281	179.463	180.196	180.451	180.374	1.2

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-o-o-

เอกสารไม่ควบคุม

เอกสารไม่ควบคุม

FOSS

Customer Service Report

Date:	Jan 27, 2025	Report No.:	13319
Job No.:	11675	Customer:	UAE
Instrument:	KT200	Address:	Bangkok
		Serial:	91790524
Start	Travel To Customer (Hrs)	Labour (Hrs)	Travel From Customer (Hrs)
Finish	08:00 10:00	10:00 13:00	1 3

Application	Special	Standard
Distributor	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Digital Service	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Internal	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Investigate	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
	Courtesy Visit	Installation
	PMA Onboarding	Quote
	Warranty	Repair
	Sales Support	Remote
		Health Check Visit

PMA Type	Smartcare	Smartcare Pro	Fosscore
	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
	Smartcare Advance	Fosscore Pro	N/A

Details of Work / Test				
- PM -				
+ Visual Check +				
- No leak				
- have damage on heater & main switch				
- ok				
- Rot of				
+ below heater 10. main switch information via				
- ok				
+ below PM kit #1 set				
- ok				
+ Function Check				
- Power on/off				
- Alkal				
- Steam				
- Condenser				
- Water pump				
- ok				
ok				
Instrument Ready for Use	OK	in	Not OK*	x

Part No.	Batch	Description	Qty
10069965	11.08.2024	Foss PM kit KT200 8.5ml Analyser/210g	1
10003912	28.03.2024	Test in element stream	1
15630111	19.10.2022	Sm. fda R595kit + 2 Pa	1

I confirm this report is accurate and complete	
Signed FOSS	Signed Customer
Name	Name
Email	Customer Contact:
*Remark: เอกสารไม่ควบคุม	

FOSS

Customer Service Report

Date:	24 February 2025	Report No.:	13854
Job No.:	11735	Customer:	UAE
Instrument:	KT8100	Address:	Bangkok
		Serial:	91887052
Start	Travel To Customer (Hrs)	Labour (Hrs)	Travel From Customer (Hrs)
Finish	07:00 09:00	09:00-12:00 13:00-14:00	1 0hrs

Application	Special	Standard
Distributor	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Digital Service	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Internal	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Investigate	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
	Courtesy Visit	Installation
	PMA Onboarding	Quote
	Warranty	Repair
	Sales Support	Remote
		Health Check Visit

PMA Type	Smartcare	Smartcare Pro	Fosscore
	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
	Smartcare Advance	Fosscore Pro	N/A

Details of Work / Test				
#PM KT200 12 Mo				
- test before 0m				
- cleaning KT200, 36 Mo replace				
- flushing Alkali RMP				
- test operation				
- Distillation 80 - 80 ml				
- Distillation 6 min 150-170ml				
- Alkali 50-80ml				
- All pass				

Part No.	Batch	Description	Qty
60031810	08.01.2024	Foss PM kit KT8100 8.100 3.6 Mo	1

I confirm this report is accurate and complete	
Signed FOSS	Signed Customer
Name	Name
Email	Customer Contact:
*Remark: เอกสารไม่ควบคุม	

Please scan QR code



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert.No.: 25CH263
Page.: 1 of 3

Equipment : pH Meter
Manufacturer : Horiba
Model : LAQUA-PH210
Serial No. : HA1M0043
ID No. : UAE.EFM.013/2565(EFM.pH.03/65)
Condition As-Received: Used Item
Received Date : 25 February 2025
Calibration Date : 26 to 28 February 2025
Reference : 2502-0783WSC-3
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phraekhanong, Bangkok 10260
Ambient Temperature : (25 ± 2.5) °C
Relative Humidity : (50 ± 15) %
Calibration Procedure : In-house method :
- CP-CH5 by direct measurement with DC voltage
standard and direct measurement with
certified reference material (CRM)
- CP-CH8 by comparison with temperature standard
Calibrated by : Warakorn Lernagatrakul
Approved by :
() Chakrit Waewwanjua
() Ponpan Paipim
(✓) Saitip Meangmai
Issue Date : 28 February 2025

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3: Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Cert.No.: 25CH263
Page.: 2 of 3

Condition of this calibration result

Instrument	Serial No.	ID No.	Cert. No.	Due Date
1) Document Process Calibrator	54030049	130RC116	24E2759	25 Aug 2025
2) Ref. Standard Thermometer	4982054	110RC044	24I757	14 July 2025

- This Certification is traceable to SI Through Technology Promotion Association (Thailand - Japan)

2. Certified Reference Materials : The measurement results are traceable to SI through Hach Lenge GmbH Ltd., Deutsche Akkreditierungsstelle, Accredited No.D-RM-15184-01-00
: The measurement results are traceable to SI through CPA chem Ltd., ANSI-ASQ National Accreditation Board, Accredited No. AR-1835

Buffer Solution	Manufacturer	Lot No.	Exp. date
pH 4.007	CPA chem	1066665	18 Jan 2027
pH 6.999	Hach Lenge GmbH	C03220	29 Oct 2026
pH 10.010	CPA chem	1066669	18 Jan 2026

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration Results

Function : mV Measurement

Performing standard curve by Document Process Calibrator at pH (4.7)(7.10)

Unit Under Calibration	Nominal Value	Standard Voltage Input	Actual Reading	Uncertainty of Measurement	Coverage factor
	pH	mV	mV	(mV)	k
pH Meter	4.00	177.48	177.4	0.058	2.00
S/N: HA1M0043	7.00	0.00	-0.1	0.058	2.00
	7.00	0.00	-0.1	0.058	2.00
	10.00	-177.48	-177.6	0.058	2.00

เอกสารไม่ควบคุม



Calibration Results

Function : pH Measurement

Performing three buffers standard curve by using buffer nominal pH (4,7)(7,10)

Unit Under Calibration	Standard pH Buffer Solution	Actual pH Reading	Actual mV Reading (mV)	Uncertainty of pH Measurement (\pm)	Coverage factor k
pH Electrode S/N.: 992H0385	4.007	4.01	147.9	0.0085	2.05
	6.999	7.00	-24.3	0.0092	2.00
	6.999	7.01	-24.4	0.0085	2.00
	10.010	10.01	-197.8	0.0092	2.00

Function : Temperature Measurement

(*) Without adjustment

This equipment was connected with Temperature Probe;

- Model : 9652
- Serial No. : 992H0385

Dimension of probe

- Length : 110 mm.
- Diameter : 16 mm.
- Immersion Depth : 80 mm.

Calibration Point (°C)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty of measurement (\pm °C)	Coverage factor k
15.0	15.003	15.0	-0.003	0.13	2.00
30.0	30.004	30.0	-0.004	0.13	2.00
45.0	45.002	45.0	-0.002	0.13	2.00

Remark - UUC* = Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k , providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่ควบคุม

ภาคผนวก ซ

หนังสือขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์

ที่ อก ๐๓๑๐(๑)/ ๘๓ ๐ ๒



กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๐๘ ตุลาคม ๒๕๖๘

เรื่อง ยกเลิกบุคลากรของห้องปฏิบัติการวิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท ยูไนเต็ด แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๒๕ กันยายน ๒๕๖๘

ตามคำขอที่อ้างถึง บริษัท ยูไนเต็ด แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด
ห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ๖-๑๔๕ สถานที่ตั้งเลขที่ ๓ ซอยอุดมสุข ๔๑ ถนนสุขุมวิท
แขวงบางจาก เขตพระโขนง กรุงเทพมหานคร ขอยกเลิกบุคลากร ความละเอียดแจ้งแล้ว นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้ออกเลิกเจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกชน
จำนวน ๓ ราย ได้แก่

- | | |
|--------------------------|----------------------------|
| ๑) นางสาวอารียา ทารามย์ | ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๐๖๗ |
| ๒) นางสาวศิริเพชร ทองขาว | ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๑๑๗ |
| ๓) นางสาวปวีณา แดนชนบท | ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๑๖๔ |

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

(นางสาวปัทมวรรณ คุณประเสริฐ)
ผู้อำนวยการกองวิจัยและเฝ้าระวังมลพิษโรงงาน
ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม



ดำเนินการถูกต้อง

กองวิจัยและเฝ้าระวังมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๕๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕

โทรสาร ๐ ๒๕๓๐ ๖๓๑๒ ต่อ ๒๑๕๔

ไปรษณีย์อิเล็กทรอนิกส์ saraban@diw.mail.go.th



"อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว"



ที่ อก ๐๓๑๐(๑)/ ๕๖ ๙ ๑

กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๐๗ กรกฎาคม ๒๕๖๘

เรื่อง เปลี่ยนแปลงบุคลากรและสารมลพิษที่วิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท ยูไนเต็ด แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๙ พฤษภาคม ๒๕๖๘

สิ่งที่ส่งมาด้วย เอกสารแนบท้ายหนังสือเปลี่ยนแปลงบุคลากรและสารมลพิษที่วิเคราะห์
บริษัท ยูไนเต็ด แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด จำนวน ๒ แผ่น

ตามหนังสือที่อ้างถึง บริษัท ยูไนเต็ด แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด
ห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ๖-๑๔๕ สถานที่ตั้งเลขที่ ๓ ซอยอุดมสุข ๔๑ ถนนสุขุมวิท
แขวงบางจาก เขตพระโขนง กรุงเทพมหานคร ขอเปลี่ยนแปลงบุคลากรและสารมลพิษที่วิเคราะห์ที่
ความละเอียดแจ้งแล้ว นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว มีความเห็นดังนี้

๑. ให้ออกเลิกเจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกชน จำนวน ๓ ราย

- | | |
|------------------------------|----------------------------|
| ๑) นายสุสันต์ พันสิทธิ์ | ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๐๐๑ |
| ๒) นางสาวสุภัทสร่า เอียนเงิน | ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๑๔๔ |
| ๓) นางสาวชานันดา กิมาคม | ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๑๘๖ |

๒. ให้เพิ่มผู้ควบคุมห้องปฏิบัติการวิเคราะห์เอกชน จำนวน ๑ ราย

นายสุสันต์ พันสิทธิ์

ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๐๔๗

๓. ให้เพิ่มข้อข่ายสารมลพิษที่วิเคราะห์ในน้ำใต้ดิน อากาศเสีย และดิน ตามสิ่งที่ส่งมาด้วย
อนึ่ง หนังสือฉบับนี้จะมีผลตั้งแต่วันที่ออกให้ขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
ในวันที่ ๒ กุมภาพันธ์ ๒๕๖๘

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

(นายประสม คำวังพงษ์)
ผู้อำนวยการกองวิจัยและเฝ้าระวังมลพิษโรงงาน
ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม

กองวิจัยและเฝ้าระวังมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๕๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕

โทรสาร ๐ ๒๕๓๐ ๖๓๑๒ ต่อ ๒๑๕๔

ไปรษณีย์อิเล็กทรอนิกส์ saraban@diw.mail.go.th



"อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว"



ที่ อก ๐๓๑๐(๑)/ ๖๙ ๗ ๓



กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๒๕ สิงหาคม ๒๕๖๘

เรื่อง ยกเลิกบุคลากรของห้องปฏิบัติการวิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท ยูไนเต็ด แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๗ สิงหาคม ๒๕๖๘

ตามคำขอที่อ้างถึง บริษัท ยูไนเต็ด แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด
ห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ๖-๑๔๕ สถานที่ตั้งเลขที่ ๓ ซอยอุดมสุข ๔๑ ถนนสุขุมวิท
แขวงบางจาก เขตพระโขนง กรุงเทพมหานคร ขอยกเลิกบุคลากร ความละเอียดแจ้งแล้ว นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้ออกเลิกเจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกชน
จำนวน ๕ ราย ได้แก่

- | | |
|---------------------------|----------------------------|
| ๑) นายคณิติน พงษ์ศรีพานุร | ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๑๐๔ |
| ๒) นายธีรวัฒน์ ธรรมสุวรรณ | ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๑๓๗ |
| ๓) นายอาทิตย์ ตากา | ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๑๖๒ |
| ๔) นางสาวกมลชนก ปูนคำ | ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๑๓๒ |
| ๕) นายวีระพงษ์ แสงทำนัง | ทะเบียนเลขที่ ๖-๑๔๕-๖-๐๑๘๘ |

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

(นายศิระ จันทร์ลิล)
นักวิทยาศาสตร์เชี่ยวชาญ วิชาการจัดการแผน
ผู้อำนวยการกองวิจัยและเฝ้าระวังมลพิษโรงงาน
ปฏิบัติราชการแทนอธิบดีกรมโรงงานอุตสาหกรรม

กองวิจัยและเฝ้าระวังมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๕๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕

โทรสาร ๐ ๒๕๓๐ ๖๓๑๒ ต่อ ๒๑๕๔

ไปรษณีย์อิเล็กทรอนิกส์ saraban@diw.mail.go.th



ดำเนินการถูกต้อง



"อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว"



เอกสารแนบท้ายหนังสือเปลี่ยนแปลงบุคลากรและสารมลพิษที่วิเคราะห์
บริษัท ยูไนเต็ด แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด เลขทะเบียน ๖-๑๔๕
ที่ อก ๐๓๑๐(๑)/ ๕๖ ๙ ๑ ลงวันที่ ๐๗ กรกฎาคม ๒๕๖๘

ขอขยายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๑๑ รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Aluminum	Digestion, Inductively Coupled Plasma Method ^[1]
2	Copper	1) Digestion, Flame Atomic Absorption Spectrometric Method ^[1] 2) Digestion, Inductively Coupled Plasma Method ^[1]
3	Iron	1) Digestion, Flame Atomic Absorption Spectrometric Method ^[1] 2) Digestion, Inductively Coupled Plasma Method ^[1]
4	Molybdenum	Digestion, Inductively Coupled Plasma Method ^[1]

อากาศเสีย (ปล่อยระบาย) จำนวน 1 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Oxides of Nitrogen	Absorption Sampling, Ion Chromatographic Method ^[2]

ดิน จำนวน 6 รายการ

ลำดับที่	สารมลพิษ	วิธีวิเคราะห์
1	Aluminum	Digestion, Inductively Coupled Plasma Method ^[3,5]
2	Copper	1) Digestion, Flame Atomic Absorption Spectrometric Method ^[3,6] 2) Digestion, Inductively Coupled Plasma Method ^[3,5]
3	Iron	1) Digestion, Flame Atomic Absorption Spectrometric Method ^[3,6] 2) Digestion, Inductively Coupled Plasma Method ^[3,5]
4	Molybdenum	Digestion, Inductively Coupled Plasma Method ^[3,6]
5	pH	Electrometric Method ^[8]
6	TPH (C ₅ -C ₈)	Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^[4,7]

เอกสารอ้างอิง

1. APHA, AWWA, WEF. Standard Methods for the Examination of Water and Wastewater. 24th ed. Washington, DC: APHA, 2023.
2. United States Environmental Protection Agency. Standards of Performance for New Stationary Sources. 40 CFR 60. Appendix A, 2023.

3. United States...

3. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Acid Digestion of Sediments, Sludges, and Soils. SW-846 Method 3050B, 1996.
4. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Volatile Organic Compounds in Various Sample Matrices Using Equilibrium Headspace Analysis. SW-846 Method 5021A, 2014.
5. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Inductively Coupled Plasma-Optical Emission Spectrometry. SW-846 Method 6010D, 2018.
6. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Flame Atomic Absorption Spectrophotometry. SW-846 Method 7000B, 2007.
7. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Volatile Organic Compounds by Gas Chromatography/Mass Spectrometry. SW-846 Method 8260D, 2018.
8. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Soil and Waste pH. SW-846 Method 9045D, 2004.

กาน์



ที่ อก ๐๓๑๐(๑)/ ๑ ๘ ๘ ๑

กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๒๔ กุมภาพันธ์ ๒๕๖๘

เรื่อง ยกเลิกบุคลากรของห้องปฏิบัติการวิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท ยูไนเต็ด แอนาไลติกส์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๑๘ กุมภาพันธ์ ๒๕๖๘

ตามคำขอที่อ้างถึง บริษัท ยูไนเต็ด แอนาไลติกส์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด
ห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ๖-๑๔๕๕ สถานที่ตั้งเลขที่ ๓ ซอยอุดมสุข ๔๑ ถนนสุขุมวิท
แขวงบางจาก เขตพระโขนง กรุงเทพมหานคร ขอยกเลิกบุคลากร ความละเอียดแจ้งแล้ว นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้ยกเลิกเจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกชน
จำนวน ๓ ราย ได้แก่

- | | |
|---------------------------|-----------------------------|
| ๑) นายอภิสิทธิ์ ศรีคงแก้ว | ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๕๘ |
| ๒) นางสาวนันธิดา พรหมฉาย | ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๕๙ |
| ๓) นายภูวดล เปิงมา | ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๕๘ |

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

รังสรวง
(นายธีรทัศน์ อิศรางกูร ณ อยุธยา)
รองอธิบดี ปฏิบัติราชการแทน
อธิบดีกรมโรงงานอุตสาหกรรม

กองวิจัยและเตือนภัยมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕

โทรสาร ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๔๔

ไปรษณีย์อิเล็กทรอนิกส์ sarabang@diw.mail.go.th



UAE
UNITED ANALYST AND ENGINEERING
CONSULTANT COMPANY LIMITED
สำนักงานลูกค้า

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ กองวิจัยและเตือนภัยมลพิษโรงงาน กรมโรงงานอุตสาหกรรม โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕



"อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว"



ที่ อก ๐๓๑๐(๑)/ ๑ ๐ ๘ ๘

กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๖ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๔๐๐

๐๗ กุมภาพันธ์ ๒๕๖๘

เรื่อง ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

เรียน กรรมการผู้จัดการ บริษัท ยูไนเต็ด แอนาไลติกส์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ต่ออายุ/เปลี่ยนแปลงบุคลากร และชนิดสารมลพิษของห้องปฏิบัติการวิเคราะห์เอกชน
ลงวันที่ ๓ ธันวาคม ๒๕๖๗

- สิ่งที่ส่งมาด้วย ๑. รายชื่อผู้ควบคุมห้องปฏิบัติการวิเคราะห์เอกชน จำนวน ๔๐ ราย
๒. รายชื่อเจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกชน จำนวน ๑๔๑ ราย
๓. ขอบข่ายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม

ตามคำขอที่อ้างถึง บริษัท ยูไนเต็ด แอนาไลติกส์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด ขอต่ออายุ
หนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน เลขทะเบียน ๖-๑๔๕๕ สถานที่ตั้งเลขที่ ๓ ซอยอุดมสุข ๔๑
ถนนสุขุมวิท แขวงบางจาก เขตพระโขนง กรุงเทพมหานคร ต่อกรมโรงงานอุตสาหกรรม นั้น

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้บริษัท ยูไนเต็ด แอนาไลติกส์ แอนด์ เอ็นจิเนียริง
คอนซัลแตนท์ จำกัด ต่ออายุหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน โดยมีองค์ประกอบดังนี้

- ก. ผู้ควบคุมห้องปฏิบัติการวิเคราะห์เอกชน จำนวน ๔๐ ราย ตามสิ่งที่ส่งมาด้วย ๑
ข. เจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกชน จำนวน ๑๔๑ ราย ตามสิ่งที่ส่งมาด้วย ๒
ค. ขอบข่ายสารมลพิษที่ได้รับขึ้นทะเบียนให้วิเคราะห์ในน้ำ/น้ำเสีย น้ำใต้ดิน อากาศเสีย
สิ่งปฏิกูลหรือวัสดุที่ไม่ใช้แล้ว และดิน ตามสิ่งที่ส่งมาด้วย ๓

หนังสือฉบับนี้จะหมดอายุในวันที่ ๒ กุมภาพันธ์ ๒๕๗๒ หากประสงค์จะต่ออายุหนังสือ
รับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน ให้ยื่นคำขอต่ออายุพร้อมเอกสารประกอบคำขอต่อกรมโรงงาน
อุตสาหกรรมภายใน ๖๐ วัน ก่อนวันสิ้นสุดอายุของหนังสือรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

รังสรวง
(นายธีรทัศน์ อิศรางกูร ณ อยุธยา)
รองอธิบดี ปฏิบัติราชการแทน
อธิบดีกรมโรงงานอุตสาหกรรม



UAE
UNITED ANALYST AND ENGINEERING
CONSULTANT COMPANY LIMITED
สำนักงานลูกค้า

กองวิจัยและเตือนภัยมลพิษโรงงาน

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและทะเบียนห้องปฏิบัติการ

โทร. ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๓-๕

โทรสาร ๐ ๒๔๓๐ ๖๓๑๒ ต่อ ๒๑๐๔๔

ไปรษณีย์อิเล็กทรอนิกส์ sarabang@diw.mail.go.th



"อุตสาหกรรมก้าวไกล ประเทศไทยก้าวหน้า ร่วมกันพัฒนา อุตสาหกรรมสีเขียว"



เอกสารแนบท้ายหนังสือต่ออายุรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน

บริษัท ยูไนเต็ด แอนาไลติกส์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด เลขทะเบียน ๖-๑๔๕๕

ที่ อก ๐๓๑๐(๑)/ ๑ ๐ ๘ ๘

ลงวันที่ ๐๗ กุมภาพันธ์ ๒๕๖๘

สิ่งที่ส่งมาด้วย ๑

ก. ผู้ควบคุมห้องปฏิบัติการวิเคราะห์เอกชน จำนวน ๔๐ ราย

๑) นางสาวกฤษกรณ กัทธิกรกุล

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๐๑

๒) นายณรงค์ อิมพาลี

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๐๒

๓) นางสาวนันธิดา บุญไชย

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๐๓

๔) นางปิยะพัชร สุทธิมนัสวงษ์

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๐๔

๕) นางสาวเบญจวรรณ วิริยชัย

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๐๖

๖) นายพรรัตน์ วงศ์บุรุษชัย

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๐๗

๗) นางสาวฉวีวรรณ บุญลา

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๐๘

๘) นายสุวิทย์ จอดนอก

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๐๙

๙) นางสาวโชติกา สมบูรณ์

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๑๐

๑๐) นางสาวบุษกร เลิศกาญจนา

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๑๑

๑๑) นางสาววิไลลักษณ์ ศรีสุข

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๑๒

๑๒) นายศิลา บรรจงใจรักษ์

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๑๓

๑๓) นายปฏิกรณ์ คณะนา

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๑๔

๑๔) นายธีรวัฒน์ ขมิ้ม

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๑๕

๑๕) นางสาวศิริพร ศรีประดิษฐ์

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๑๖

๑๖) นางสาวศิริพร ศรีประดิษฐ์

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๑๗

๑๗) นางสาวศิริพร ศรีประดิษฐ์

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๑๘

๑๘) นายภูษงค์ พาณิชย์เลิศอาพา

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๑๙

๑๙) นายณัฐวัฒน์ แดงสวัสดิ์

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๒๐

๒๐) นายเอกรัตน์ ปะกะนิมิต

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๒๑

๒๑) นางสาวนิศารัตน์ ศรีสกุลสิงห์โชค

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๒๒

๒๒) นางสาวเจตจิรพร ทวีระอด

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๒๓

๒๓) นางสาวสุวรรณ คงทอง

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๒๔

๒๔) นางสาววรรณ พัดสองชั้น

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๒๕

๒๕) นายวิรัชพงศ์ โภกแก้ว

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๒๖

๒๖) นายวัชรพงษ์ เทพดนตรี

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๒๗

๒๗) นายอนุชา สยดี

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๒๘

๒๘) นายกรวิทย์ เชื้อศิริกุล

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๒๙

๒๙) นายสุธีระ อรุณจันทร์

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๓๐

๓๐) นางสาวทัศนีย์ อ่อนคำ

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๓๑

๓๑) นางพริ้มพรรณ กองสิน

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๓๒

๓๒) นายศุภณัฐ คุณอนาญจ

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๓๓

๓๓) นางสาวศิริพร เหมือนเ

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๓๔

๓๔) นางสาวณิลล์ ขำนิล

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๓๕

๓๕) นางสาวพริมา ธีระจินดา

ทะเบียนเลขที่ ๖-๑๔๕๕-๖-๐๐๓๖

กาน์

๓๖) นายนาเคนทร์...

๓๖) นายนาเคนทร์ พันธุ์ชาติกุล
๓๗) นายกานต์พงศ์ บุญพวง
๓๘) นางสุธรรมมา แก้วชื่อนอก
๓๙) นางสาวสริน โขยเชษฐ์พิพัฒกุล
๔๐) นางมานิดา แยมโย

ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๔๐
ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๔๑
ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๔๒
ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๔๓
ทะเบียนเลขที่ ๖-๑๔๕-ค-๐๐๔๔

๑๖๖



สำนักงานที่ดิน

เอกสารแนบท้ายหนังสือคำอุทธรณ์ขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกสาร
บริษัท ยูไนเต็ด แอแนลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด เลขทะเบียน ๖-๑๔๕-
ที่ อก ๐๓๑๐(๑) / ๑๐๘๔ ลงวันที่ ๐๗ กุมภาพันธ์ ๒๕๖๘

ข. เจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์เอกสาร จำนวน ๑๔๑ ราย

๑) นายสุสันต์ พันธุ์สิงห์
๒) นายพัณณัฐ เจริญผล
๓) นางสาววิไลลักษณ์ เกาโลง
๔) นายสมชาติ อุทุมรัตน์
๕) นางสาวปรมาภรณ์ ทองแก้ว
๖) นางสาวกัญญา สมพงษ์
๗) นางสาววรรณิ์ สายบุญเรือน
๘) นายกฤษณพงษ์ นามทิพย์
๙) นางสาวอารณีย์ อ่อนคง
๑๐) นายกิตติศักดิ์ ทรงจำรัส
๑๑) นางสาวอักษรินทร์ บุญคง
๑๒) นางสาวพรพิมล แวนทอง
๑๓) นายอภิรักษ์ ทั่วทั้ง
๑๔) นายมานิตย์ ปานโชติ
๑๕) นายทศพร ธนะพิรุฬห์
๑๖) นางสาวกัญญาณี โยธา
๑๗) นางสาวเกสรี สุทธิ
๑๘) นางสาวชนอนันท์ อภิพัทธ์ปภา
๑๙) นายศิริพัชร จงผดุงเกียรติ
๒๐) นางสาวสุภาวดี อินยาศิริ
๒๑) นายพงศ์เทพ เหล่าจรร
๒๒) นายขวัญชัย พันธุ์
๒๓) นางสาวพัชจิรา คดีพิศาล
๒๔) นางสาวเมกกา เลิศจำนงค์
๒๕) นายพีระพัฒน์ บุญญศิริศิลป์
๒๖) นายชัชวาลย์ เลื่อนทอง
๒๗) นายณัฐกร ธรรมรัตน์
๒๘) นายกันนิกร ระโส
๒๙) นายปริญญา กลมเกลียว
๓๐) นายธีรวัฒน์ มาตราโพธิ์ศรี
๓๑) นายบุญฤทธิ์ ก้อนสิน
๓๒) นายพรชภูมิ โกวสกุล
๓๓) นายชิตะ แสงจันทร์
๓๔) ว่าที่ร้อยตรีณัฐพงศ์ เมืองชัย
๓๕) นายอนันท์ เลิศประเสริฐ

ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๐๑
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๐๒
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๐๓
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๐๔
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๐๕
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๐๖
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๐๗
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๐๘
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๐๙
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๑๐
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๑๑
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๑๒
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๑๓
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๑๔
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๑๕
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๑๖
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๑๗
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๑๘
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๑๙
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๒๐
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๒๑
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๒๒
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๒๓
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๒๔
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๒๕
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๒๖
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๒๗
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๒๘
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๒๙
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๓๐
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๓๑
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๓๒
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๓๓
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๓๔
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๓๕
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๓๖
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๓๗
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๓๘
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๓๙
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๐
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๑
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๒
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๓
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๔
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๕
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๖
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๗
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๘
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๙
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๐

๑๖๖

๓๖) นางสาวนิภาพร...

๓๖) นางสาวนิภาพร จันทร์เขตต์
๓๗) นายรณภพ ภูตระกูลพัฒนา
๓๘) นายสมพงษ์ สกุลไทย
๓๙) นายสุวิทย์ นิธิจิตชูวงศ์
๔๐) นายอัยยกุล วนศิริ
๔๑) นายเอกภูมิ เสนอใจ
๔๒) นายสุสันต์ บุญเลี้ยง
๔๓) นายธนเดช หวานเสนาะ
๔๔) นายอภิสิทธิ์ ศรีคงแก้ว
๔๕) ว่าที่ร้อยตรีอุทัย แก้วรามภู
๔๖) นางสาวนารินทร์ สานนท์
๔๗) นายศุภกร รินวงศ์
๔๘) นางสาวจินตสุภา เป็เลียนศรี
๔๙) นางสาวเนตรนภา กมลบุรินทร์
๕๐) นางสาวอารียา ทรากรมย์
๕๑) นายจิรวัฒน์ สุขเกษม
๕๒) นายกิตติพงษ์ สอนชัยภูมิ
๕๓) นายจุมพล สอนเพชร
๕๔) นางสาวพัชรภรณ์ แสงฟ้า
๕๕) นายรัตนชัย เหล่ามา
๕๖) นายอิทธิพงษ์ ศรีวิเศษ
๕๗) นางสาวกรรณิการ์ สำสันท
๕๘) นางสาวณัฐชา พรหมศิริ
๕๙) นายณัฐสิทธิ์ ศรีพิมพ์
๖๐) นางสาวลักขณา จันทร์สุข
๖๑) นายกิตติศักดิ์ นุ่มมิม
๖๒) นายวรพงษ์ นพจันทร์
๖๓) นางสาวชนาภา มาคะมาตร
๖๔) นายณัฐชัย พรหมอารักษ์
๖๕) นายณรินทร์ พายแก้ว
๖๖) นายปรีชาพล โสภา
๖๗) นายวิรัตน์ แสนงาม
๖๘) นายอาทิตย์ อุดมผล
๖๙) นายอิทธิเดช ใจบุญ
๗๐) นายชิตานันท์ พงษ์จิตรานุสร
๗๑) นายเสกสรรค์ เอกกลิ่นบัว
๗๒) นางสาวนันทา แหวนในเมือง
๗๓) นางสาวพิมลวรรณ สิมมา

ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๗
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๘
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๔๙
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๐
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๑
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๒
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๓
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๔
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๕
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๖
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๗
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๘
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๕๙
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๖๐
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๖๑
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๖๒
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๖๓
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๖๔
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๖๕
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๖๖
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๖๗
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๖๘
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๖๙
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๗๐
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๗๑
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๗๒
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๗๓
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๗๔
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๗๕
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๗๖
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๗๗
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๗๘
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๗๙
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๘๐
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๘๑
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๘๒
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๘๓
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๘๔
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๘๕
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๘๖
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๘๗
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๘๘
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๘๙
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๙๐
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๙๑
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๙๒
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๙๓
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๙๔
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๙๕
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๙๖
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๙๗
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๙๘
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๐๙๙
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๐๐
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๐๑
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๐๒
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๐๓
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๐๔
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๐๕
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๐๖
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๐๗
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๐๘
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๐๙
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๑๐

๑๖๖

๓๗) นายณัฐวัฒน์...

๓๕) นายณัฐวัฒน์ วงศ์คำ
๓๖) นายประพันธ์ยุทธ เลื่อนนาง
๓๗) นางสาวกัญญา ลำตัน
๓๘) นางสาวนิภาพร ชื่นนุกัม
๓๙) นางสาวบุญญา มอญคุณ
๔๐) นายอมรพร อรรถลักษณ์
๔๑) นางสาวศรีเพียร ทองขาว
๔๒) นางสาวนิภาพร ศุภาดิโกสร
๔๓) นางสาววิมลวรรณ คำตัน
๔๔) นายคุณานันท์ ฤทธากานานนท์
๔๕) นายชาญณรงค์ อัครชัย
๔๖) นางสาวจิตติมา ศรีวรรณ
๔๗) นายสุจิตต์ ปิ่นเงิน
๔๘) นายเชษฐา ชัยศิริ
๔๙) นายชัช หนองบัว
๕๐) นายชัย บัวสรวง
๕๑) นางสาวอรุณา ประสานศรี
๕๒) นายพดล เบียมเนียม
๕๓) นายศุภกร สอนศรี
๕๔) นายคุณพล ศิลานนท์
๕๕) นายโชคชัย พุ่มไสว
๕๖) นายธีรวัฒน์ ธรรมสุวรรณ
๕๗) นายณัฐพงศ์ ชะขุนทด
๕๘) นางสาวณัฐกานดา พลนิกรกิจ
๕๙) นางสาวพนพร ทองบุญ
๖๐) นางสาวพรชิตา ขจรนดิษฐ์
๖๑) นางสาวณัฐกานดา รอดทอง
๖๒) นางสาวณัฏฐา แสงสว่าง
๖๓) นายกริณี สีใส
๖๔) นายศุภพร คงศรี
๖๕) นางสาวสุกัลยา เอี่ยมเงิน
๖๖) นางสาวพรพรรณ ฟ้าอ่อน
๖๗) นายอนันต์ นาค
๖๘) นางสาวพรพิมล ประชาพัชร์
๖๙) นายธีรภัทร บุญยง
๗๐) นางสาวณัฐชา แก้วภาพ
๗๑) นายสิทธิพล พร้อมพงษ์บุญ
๗๒) นางสาวนันทา กลิ่นหนู

ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๑๑
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๑๒
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๑๓
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๑๔
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๑๕
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๑๖
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๑๗
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๑๘
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๑๙
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๒๐
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๒๑
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๒๒
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๒๓
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๒๔
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๒๕
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๒๖
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๒๗
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๒๘
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๒๙
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๓๐
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๓๑
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๓๒
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๓๓
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๓๔
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๓๕
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๓๖
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๓๗
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๓๘
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๓๙
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๔๐
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๔๑
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๔๒
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๔๓
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๔๔
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๔๕
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๔๖
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๔๗
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๔๘
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๔๙
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๕๐
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๕๑
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๕๒
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๕๓
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๕๔
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๕๕
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๕๖
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๕๗
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๕๘
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๕๙
ทะเบียนเลขที่ ๖-๑๔๕-จ-๐๑๖๐

๑๖๖

๑๓๖) นางสาวปิยา...

- ๔ -

๑๓๓) นางสาวปติยา ชูเชิดเชื้อ
๑๓๔) นางสาวลัดดาวัลย์ โพธิ์พันธ์
๑๓๕) นายอาทิตย์ ตาภา
๑๓๖) นางสาวบุญยาพร บุญอมศรี
๑๓๗) นางสาวพัชรารวรรณ จันธิบุตร
๑๓๘) นางสาวนฤกร ใต้น้ำกวย
๑๓๙) นางสาวปวีณา แดนชนบท
๑๔๐) นางสาวนันทิศา พรหมกวด
๑๔๑) นางสาวกมลชนก ปูนคำ
๑๔๒) นางสาวปาริฉัตร ทองใบ
๑๔๓) นายชัยวัฒน์ จันละคร
๑๔๔) นางสาวกัญญา สิงห์แก้ว
๑๔๕) นางสาวอารีนา มณีเชีย
๑๔๖) นายสุภากรณ์ อนุรา
๑๔๗) นางสาวชานันดา กิมาคม
๑๔๘) นายอนันต์นรินทร์ ยาสี
๑๔๙) นายวีระพงษ์ แสงทำนง
๑๕๐) นางสาวปิยะนัฐชา สำเนาพงษ์
๑๕๑) นางสาวนันทิสรา ศรีสถาน
๑๕๒) นางสาวจุฑิษัณธ์ โสแทน
๑๕๓) นายธีรวัฒน์ พรหมลา
๑๕๔) นายอนันต์ชัย ปลั่งกลาง
๑๕๕) นายณภัทร เตมียบุตร
๑๕๖) นางสาวจิตาภา ฤาชา
๑๕๗) นางสาวสมณฑาทิพย์ สังข์ทอง
๑๕๘) นางสาววชิรา บาบัญ
๑๕๙) นายภูวดล เบื้องมา
๑๖๐) ว่าที่ร้อยตรีณยุทธ ประทุมเขตต์
๑๖๑) นายอนุสรณ์ พลสำโรง

ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๖๐
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๖๑
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๖๒
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๖๓
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๖๔
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๖๕
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๖๖
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๖๗
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๖๘
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๖๙
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๗๐
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๗๑
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๗๒
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๗๓
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๗๔
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๗๕
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๗๖
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๗๗
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๗๘
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๗๙
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๘๐
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๘๑
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๘๒
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๘๓
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๘๔
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๘๕
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๘๖
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๘๗
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๘๘
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๘๙
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๙๐
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๙๑
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๙๒
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๙๓
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๙๔
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๙๕
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๙๖
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๙๗
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๙๘
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๖๙๙
ทะเบียนเลขที่ ๖-๑๔๕-๖-๑๗๐๐

UAE
UNITED ANALYST AND ENGINEERING
CONSULTANT COMPANY LIMITED

ดำเนินการถูกต้อง

เอกสารแนบท้ายหนังสือต่ออายุรับขึ้นทะเบียนห้องปฏิบัติการวิเคราะห์เอกชน
บริษัท ยูนิเทค แอนาไลติกส์ แอนด์ เอ็นจิเนียริ่ง คอนซัลแตนท์ จำกัด เลขทะเบียน ๖-๑๔๕
ที่ ออก ๑๓๐(๑) ๑๐ ๘ ๘ ลงวันที่ ๐๗ กุมภาพันธ์ ๒๕๖๕

ขอขยายสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๓๕๗ รายการ

น้ำเสีย จำนวน 46 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Aldrin	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
2	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
3	Barium	Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
4	α-BHC	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
5	β-BHC	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
6	δ-BHC	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
7	γ-BHC	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
8	Biochemical Oxygen Demand	1) 5-Day BOD Test, Azide Modification Method ⁽⁴⁾ 2) 5-Day BOD Test, Membrane Electrode Method ⁽⁴⁾
9	Cadmium	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
10	Chemical Oxygen Demand	1) Closed Reflux, Titrimetric Method ⁽⁴⁾ 2) Closed Reflux, Colorimetric Method ⁽⁴⁾ 3) Open Reflux, Titrimetric Method ⁽⁴⁾
11	Chlordane	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
12	Chromium	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
13	Color	ADMI Weighted-Ordinate Spectrophotometric Method ⁽⁴⁾
14	Copper	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
15	Cyanide	1) Distillation, Colorimetric Method ⁽⁴⁾ 2) Total Cyanide after Distillation, by Flow Injection Analysis Method ⁽⁴⁾
16	o,p'-DDT	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
17	4,4'-DDD	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
18	4,4'-DDE	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
19	4,4'-DDT	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
20	Dieldrin	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
21	Endosulfan I	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
22	Endosulfan II	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
23	Endosulfan sulfate	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
24	Endrin	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾

25 Endrin aldehyde...

- ๒ -

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
25	Endrin aldehyde	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
26	Formaldehyde	Distillation, Colorimetric Method ⁽²⁾
27	Free Chlorine	1) Iodometric Method ⁽⁴⁾ 2) DPD Ferrous Titrimetric Method ⁽⁴⁾
28	Heptachlor	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
29	Heptachlor Epoxide	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
30	Hexavalent Chromium	Colorimetric Method ⁽⁴⁾
31	Lead	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
32	Manganese	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
33	Mercury	Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ⁽⁴⁾
34	Methoxychlor	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
35	Nickel	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
36	Oil & Grease	1) Liquid-Liquid, Partition-Gravimetric Method ⁽⁴⁾ 2) Soxhlet Extraction Method ⁽⁴⁾
37	pH	Electrometric Method ⁽⁴⁾
38	Phenols	1) Distillation, Chloroform Extraction Method ⁽⁴⁾ 2) Distillation, Direct Photometric Method ⁽⁴⁾
39	Selenium	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
40	Sulfide	1) Iodometric Method ⁽⁴⁾ 2) Methylene Blue Method ⁽⁴⁾
41	Temperature	Laboratory and Field Methods ⁽⁴⁾
42	Total Dissolved Solids	Dried at 180 °C ⁽⁴⁾
43	Total Kjeldahl Nitrogen	Semi-Micro-Kjeldahl Method ⁽⁴⁾
44	Total Suspended Solids	Dried from 103 to 105 °C ⁽⁴⁾
45	Trivalent Chromium	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾ 3) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾ 4) Colorimetric Method; Calculation ⁽⁴⁾
46	Zinc	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾

น้ำได้คืน...

- ๓ -

น้ำได้คืน จำนวน 126 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Acenaphthene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
2	Acetone	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
3	Aldrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
4	Anthracene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
5	Antimony	Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
6	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
7	Atrazine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
8	Barium	Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
9	Benz(a)anthracene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
10	Benzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
11	Benzo(b)fluoranthene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
12	Benzo(k)fluoranthene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
13	Benzoic acid	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

14 Benzo(a)pyrene...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
14	Benzo(a)pyrene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
15	Benzo(g,h,i)perylene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
16	Beryllium	Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
17	Bis(2-chloroethyl)ether	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
18	Bis(2-ethylhexyl)phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
19	Bromodichloromethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
20	Bromoform	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
21	Butanol	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
22	Butyl benzyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
23	Cadmium	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ⁽⁴⁾ 3) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
24	Carbazole	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
25	Carbon disulfide	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
26	Carbon tetrachloride	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
27	Chlordane	1) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
28	p-Chloroaniline	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

29 Chlorobenzene...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
29	Chlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
30	Chlorodibromomethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
31	Chloroform	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
32	2-Chlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
33	Chromium	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
34	Chromium (III)	1) Digestion, Direct Air-Acetylene Flame Method; Colorimetric Method; Calculation ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method; Colorimetric Method; Calculation ⁽⁴⁾
35	Chromium (VI)	Colorimetric Method ⁽⁴⁾
36	Chrysene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
37	Cyanide	Distillation, Colorimetric Method ⁽⁴⁾
38	2,4-D	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
39	DDD	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
40	DDE	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
41	DDT	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
42	Dibenz(a,h)anthracene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

43 Di-n-butyl phthalate...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
43	Di-n-butyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
44	1,2-Dichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
45	1,3-Dichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
46	1,4-Dichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
47	3,3'-Dichlorobenzidine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
48	1,1-Dichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
49	1,2-Dichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
50	1,1-Dichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
51	cis-1,2-Dichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
52	trans-1,2-Dichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
53	2,4-Dichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
54	1,2-Dichloropropane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
55	1,3-Dichloropropane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
56	1,3-Dichloropropene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
57	Dieldrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
58	Diethyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
59	2,4-Dimethylphenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
60	2,4-Dinitrophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

61 2,4-Dinitrotoluene...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
61	2,4-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
62	2,6-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
63	Di-n-Octyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
64	Endosulfan	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
65	Endrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
66	Ethylbenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
67	Fluoranthene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
68	Fluorene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
69	Heptachlor	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
70	Heptachlor epoxide	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
71	Hexachlorobenzene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
72	Hexachloro-1,3-butadiene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
73	n-Hexane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

74 α-HCH...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
74	α-HCH	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
75	β-HCH	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
76	γ-HCH	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
77	Hexachlorocyclopentadiene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
78	Hexachloroethane	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
79	Indeno(1,2,3-cd)pyrene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
80	Isophorone	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
81	Lead	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ⁽⁴⁾ 3) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
82	Manganese	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
83	Mercury	Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ⁽⁴⁾
84	Methanol	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
85	Methoxychlor	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
86	Methyl bromide	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

87 Methylene chloride...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
87	Methylene chloride	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
88	2-Methylphenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
89	2-Methylnaphthalene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
90	Methyl tert-butyl ether	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
91	Naphthalene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
92	Nickel	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
93	Nitrobenzene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
94	N-Nitrosodiphenylamine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
95	N-Nitrosodi-n-propylamine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
96	Polychlorinated Biphenyls - PCB 1016 - PCB 1221 - PCB 1232 - PCB-1242 - PCB-1248 - PCB-1254 - PCB-1260	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
97	Pentachlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
98	pH	Electrometric Method ⁽⁴⁾
99	Phenanthrene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

100 Phenol...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
100	Phenol	1) Distillation, Chloroform Extraction Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
101	Pyrene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
102	Selenium	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
103	Silver	Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
104	Styrene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
105	1,1,2,2-Tetrachloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
106	Tetrachloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
107	Toluene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
108	Toxaphene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
109	TPH (C ₅ - C ₈)	1) Purge and Trap, Gas Chromatographic Method ^(12,22) 2) Purge and Trap, Gas Chromatographic/Mass spectrometric Method ^(12,27)
110	TPH (C ₈ - C ₁₆)	Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(9,22)
111	TPH (C ₁₆ - C ₃₅)	Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(9,22)
112	1,2,4-Trichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
113	1,1,1-Trichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
114	1,1,2-Trichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
115	Trichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

116 2,4,5-Trichlorophenol...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
116	2,4,5-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
117	2,4,6-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
118	1,3,5-Trimethylbenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
119	Vanadium	Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
120	Vinyl acetate	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
121	Vinyl chloride	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
122	m-Xylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
123	o-Xylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
124	p-Xylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
125	Xylene (Total)	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
126	Zinc	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾

อากาศเสีย (ปล่องระบาย) จำนวน 25 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Antimony	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽⁵⁾
2	Arsenic	1) Isokinetic Sampling, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ⁽⁵⁾ 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽⁵⁾
3	Cadmium	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ⁽⁵⁾ 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽⁵⁾
4	Carbon Monoxide	Instrumental Analyzer Method ⁽⁵⁾
5	Chlorine	Isokinetic Sampling, Ion Chromatographic Method ⁽⁵⁾
6	Chromium	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ⁽⁵⁾

Chromium (ต่อ)...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
6	Chromium (ห่อ)	2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
7	Cobalt	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
8	Copper	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
9	Cresol	Absorption Sampling, Gas Chromatographic Method ^[5]
10	Dioxins/Furans	Isokinetic Sampling ^[5]
11	Hydrogen Chloride	Isokinetic Sampling, Ion Chromatographic Method ^[5]
12	Hydrogen Fluoride	Isokinetic Sampling, Ion Chromatographic Method ^[5]
13	Hydrogen Sulfide	Absorption Sampling, Iodometric Method ^[5]
14	Lead	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
15	Manganese	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
16	Mercury	Isokinetic Sampling, Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[5]
17	Nickel	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
18	Opacity	Ringelmann's Method ^[1]
19	Oxides of Nitrogen	1) Absorption Sampling, Phenoldisulfonic acid Method ^[5] 2) Instrumental Analyzer Method ^[5]
20	Selenium	1) Isokinetic Sampling, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[5] 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
21	Sulfur Dioxide	1) Absorption Sampling, Barium-Thorin Titrimetric Method ^[5] 2) Instrumental Analyzer Method ^[5]
22	Sulfuric Acid	Isokinetic Sampling, Barium-Thorin Titrimetric Method ^[5]

23 Total Suspended Particulate...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
23	Total Suspended Particulate	Isokinetic Sampling, Gravimetric Method ^[5]
24	Vanadium	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ^[5]
25	Xylene	1) Bag Sampling, Gas Chromatographic Method ^[5] 2) Adsorption Sampling, Gas Chromatographic Method ^[5]

สิ่งปฏิกูลหรือวัสดุที่ไม่ใช้แล้ว จำนวน 35 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Aldrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[3,9,23] 2) Ultrasonic Extraction, Gas Chromatographic Method ^[10,23]
2	Antimony	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[3,6,14] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
3	Arsenic	1) Waste Extraction, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[3,6,16] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[3,6,14] 3) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^[7,14] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]
4	Barium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[3,6,14] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
5	Beryllium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[3,6,14] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
6	Cadmium	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^[3,6,15] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[3,6,14] 3) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]
7	Chlordane	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[3,9,23] 2) Ultrasonic Extraction, Gas Chromatographic Method ^[10,23]

8 Chromium...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
8	Chromium	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^[3,6,15] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[3,6,14] 3) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]
9	Chromium (III)	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method; Waste Extraction, Colorimetric Method; Calculation ^[3,6,15,17] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method; Waste Extraction, Colorimetric Method; Calculation ^[3,6,14,17] 3) Digestion, Flame Atomic Absorption Spectrometric Method; Alkaline Digestion, Colorimetric Method; Calculation ^[7,8,15,17] 4) Digestion, Inductively Coupled Plasma Method; Alkaline Digestion, Colorimetric Method; Calculation ^[7,8,14,17]
10	Chromium (VI)	1) Waste Extraction, Colorimetric Method ^[3,17] 2) Alkaline Digestion, Colorimetric Method ^[8,17]
11	Cobalt	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[3,6,14] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
12	Copper	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^[3,6,15] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[3,6,14] 3) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]
13	2,4-D	1) Waste Extraction, Gas Chromatographic Method ^[3,9,23] 2) Ultrasonic Extraction, Gas Chromatographic Method ^[10,23]
14	DDD	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[3,9,23] 2) Ultrasonic Extraction, Gas Chromatographic Method ^[10,23]

15 DDE...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
15	DDE	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[3,9,23] 2) Ultrasonic Extraction, Gas Chromatographic Method ^[10,23]
16	DDT	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[3,9,23] 2) Ultrasonic Extraction, Gas Chromatographic Method ^[10,23]
17	Dieldrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[3,9,23] 2) Ultrasonic Extraction, Gas Chromatographic Method ^[10,23]
18	Endrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[3,9,23] 2) Ultrasonic Extraction, Gas Chromatographic Method ^[10,23]
19	Heptachlor	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[3,9,23] 2) Ultrasonic Extraction, Gas Chromatographic Method ^[10,23]
20	Lead	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^[3,6,15] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[3,6,14] 3) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]
21	Lindane	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^[3,9,23] 2) Ultrasonic Extraction, Gas Chromatographic Method ^[10,23]
22	Mercury	1) Waste Extraction, Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[3,19] 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^[3,6,14] 3) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[19] 4) Digestion, Inductively Coupled Plasma Method ^[7,14]

Mercury (ห่อ)...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
22	Mercury (ต่อ)	5) Thermal Decomposition Amalgamation and Atomic Absorption Spectrometric Method ⁽²⁰⁾
23	Methoxychlor	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
24	Molybdenum	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
25	Nickel	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^(3,6,15) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 3) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 4) Digestion, Inductively Coupled Plasma Method ^(7,14)
26	Polychlorinated Biphenyls - Aroclor 1016 - Aroclor 1221 - Aroclor 1232 - Aroclor 1242 - Aroclor 1248 - Aroclor 1254 - Aroclor 1260 - 2-Chlorobiphenyl - 2,3-Dichlorobiphenyl - 2,2',5-Trichlorobiphenyl - 2,4',5-Trichlorobiphenyl - 2,2',3,5'-Tetrachlorobiphenyl - 2,2',5,5'-Tetrachlorobiphenyl - 2,3',4,4'-Tetrachlorobiphenyl - 2,2',3,4,5'-Pentachlorobiphenyl - 2,2',4,5,5'-Pentachlorobiphenyl - 2,3,3',4,6'-Pentachlorobiphenyl - Pentachlorobiphenyl	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,24) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,24)

UAE
UNITED ANALYST AND ENGINEERING
CONSULTANT COMPANY LIMITED

ดำเนินการทดสอบ

Polychlorinated Biphenyls(ต่อ)...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
	Polychlorinated Biphenyls(ต่อ) - 2,2',3,4,4',5'-Hexachlorobiphenyl - 2,2',3,4,5,5'-Hexachlorobiphenyl - 2,2',3,5,5',6'-Hexachlorobiphenyl - 2,2',4,4',5,5'-Hexachlorobiphenyl - 2,2',3,3',4,4',5-Heptachlorobiphenyl - 2,2',3,4,4',5,5'-Heptachlorobiphenyl - 2,2',3,4,4',5,6-Heptachlorobiphenyl - 2,2',3,4,5,5',6-Heptachlorobiphenyl - 2,2',3,3',4,4',5,5',6-Nonachlorobiphenyl - Pentachlorophenol	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^(3,9,28) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28) Electrometric Method ^(31,32)
28	pH	
29	Selenium	1) Waste Extraction, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^(3,6,21) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 3) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^(7,21) 4) Digestion, Inductively Coupled Plasma Method ^(7,14)
30	Silver	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
31	Thallium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)

32 Toxaphene...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
32	Toxaphene	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
33	Trichloroethylene	1) Waste Extraction, Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(3,12,27) 2) Waste Extraction, Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(3,11,27) 3) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 4) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
34	Vanadium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
35	Zinc	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^(3,6,15) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 3) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 4) Digestion, Inductively Coupled Plasma Method ^(7,14)

ดิน จำนวน 125 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Acenaphthene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
2	Acetone	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
3	Aldrin	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
4	Anthracene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25)

Anthracene (ต่อ)...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
4	Anthracene (ต่อ)	2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
5	Antimony	Digestion, Inductively Coupled Plasma Method ^(7,14)
6	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^(7,16) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
7	Atrazine	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
8	Barium	Digestion, Inductively Coupled Plasma Method ^(7,14)
9	Benz(a)anthracene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
10	Benzene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
11	Benzo(b)fluoranthene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
12	Benzo(k)fluoranthene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
13	Benzoic acid	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
14	Benzo(a)pyrene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
15	Benzo(g,h,i)perylene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
16	Beryllium	Digestion, Inductively Coupled Plasma Method ^(7,14)

17 Bis(2-chloroethyl)ether...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
17	Bis(2-chloroethyl)ether	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
18	Bis(2-ethylhexyl)phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
19	Bromodichloromethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
20	Bromoform	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
21	Butanol	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
22	Butyl benzyl phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
23	Cadmium	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
24	Carbazole	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
25	Carbon disulfide	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
26	Carbon tetrachloride	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
27	Chlordane	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
28	p-Chloroaniline	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
29	Chlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
30	Chlorodibromomethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
31	Chloroform	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
32	2-Chlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)

33 Chromium...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
33	Chromium	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
34	Chromium (III)	1) Digestion, Flame Atomic Absorption Spectrometric Method; Alkaline Digestion, Colorimetric Method; Calculation ^(7,8,15,17) 2) Digestion, Inductively Coupled Plasma Method; Alkaline Digestion, Colorimetric Method; Calculation ^(7,8,14,17)
35	Chromium (VI)	Alkaline Digestion, Colorimetric Method ^(8,17)
36	Chrysene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
37	Cyanide	Extraction, Distillation, Colorimetric Method ^(29,30)
38	2,4-D	Ultrasonic Extraction, Gas Chromatographic Method ⁽²⁶⁾
39	DDD	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
40	DDE	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
41	DDT	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
42	Dibenz(a,h)anthracene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
43	Di-n-butyl phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
44	1,2-Dichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)

45 1,3-Dichlorobenzene...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
45	1,3-Dichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
46	1,4-Dichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
47	3,3'-Dichlorobenzidine	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
48	1,1-Dichloroethane	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
49	1,2-Dichloroethane	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
50	1,1-Dichloroethylene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
51	cis-1,2-Dichloroethylene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
52	trans-1,2-Dichloroethylene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
53	2,4-Dichlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
54	1,2-Dichloropropane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
55	1,3-Dichloropropane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
56	1,3-Dichloropropene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
57	Dieldrin	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)

58 Diethyl phthalate...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
58	Diethyl phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
59	2,4-Dimethylphenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
60	2,4-Dinitrophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
61	2,4-Dinitrotoluene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
62	2,6-Dinitrotoluene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
63	Di-n-Octyl phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
64	Endosulfan	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
65	Endrin	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
66	Ethylbenzene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
67	Fluoranthene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
68	Fluorene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
69	Heptachlor	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
70	Heptachlor epoxide	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)

Heptachlor epoxide (พี)...

ลำดับ	สารมลพิษ	วิธีการตรวจ
70	Heptachlor epoxide (คอป)	2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
71	Hexachlorobenzene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
72	Hexachloro-1,3-butadiene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
73	n-Hexane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
74	α-HCH	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
75	β-HCH	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
76	γ-HCH	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
77	Hexachlorocyclopentadiene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
78	Hexachloroethane	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
79	Indeno(1,2,3-cd)pyrene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
80	Isophorone	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
81	Lead	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
82	Manganese	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)

83 Mercury...

ลำดับตัว	สารมลพิษ	วิธีวิเคราะห์
83	Mercury	1) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^[19] 2) Thermal Decomposition Amalgamation and Atomic Absorption Spectrometric Method ^[20]
84	Methanol	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,27]
85	Methoxychlor	1) Ultrasonic Extraction, Gas Chromatographic Method ^[10,23] 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,28]
86	Methyl bromide	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,27]
87	Methylene chloride	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,27] 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^[11,27]
88	2-Methylphenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,28]
89	2-Methylnaphthalene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,28]
90	Methyl tert-butyl ether	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^[13,27]
91	Naphthalene	1) Ultrasonic Extraction, Gas Chromatographic Method ^[10,23] 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,28]
92	Nickel	1) Digestion, Flame Atomic Absorption Spectrometric Method ^[7,15] 2) Digestion, Inductively Coupled Plasma Method ^[7,14]
93	Nitrobenzene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,28]
94	N-Nitrosodiphenylamine	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,28]
95	N-Nitrosodi-n-propylamine	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^[10,28]
96	Polychlorinated Biphenyls - Aroclor 1016	1) Ultrasonic Extraction, Gas Chromatographic Method ^[10,24]

Polychlorinated Biphenyls(ต่อไป)

[illegible]

97 Pentachlorophenol...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
97	Pentachlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
98	Phenanthrene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
99	Phenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
100	Pyrene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
101	Selenium	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^(7,21) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
102	Silver	Digestion, Inductively Coupled Plasma Method ^(7,14)
103	Styrene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
104	1,1,2,2-Tetrachloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
105	Tetrachloroethylene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
106	Toluene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
107	Toxaphene	Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
108	TPH (C ₅ -C ₈)	1) Purge and Trap, Gas Chromatographic Method ^(13,27) 2) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
109	TPH (C ₉ -C ₁₆)	Ultrasonic Extraction, Gas Chromatographic Method ^(10,22)
110	TPH (C _{>16} -C ₃₅)	Ultrasonic Extraction, Gas Chromatographic Method ^(10,22)

111 1,2,4-Trichlorobenzene...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
111	1,2,4-Trichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
112	1,1,1-Trichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
113	1,1,2-Trichloroethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
114	Trichloroethylene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
115	2,4,5-Trichlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
116	2,4,6-Trichlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
117	1,3,5-Trimethylbenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
118	Vanadium	Digestion, Inductively Coupled Plasma Method ^(7,14)
119	Vinyl acetate	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
120	Vinyl chloride	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
121	m-Xylene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
122	o-Xylene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
123	p-Xylene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
124	Xylene (Total)	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)

125 Zinc...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
125	Zinc	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,13) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)

เอกสารอ้างอิง

- กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2549. เรื่อง กำหนดค่าปริมาณเคมีภัณฑ์ที่เจือปนในอากาศที่ระบายออกจากปล่องของหม้อน้ำโรงสีข้าวที่ใช้แก๊สเป็นเชื้อเพลิง. ราชกิจจานุเบกษา. 4 ธันวาคม 2549. เล่มที่ 123 ตอนพิเศษ 125ง.
- สมาคมวิศวกรรมสิ่งแวดล้อมแห่งประเทศไทย. คู่มือวิเคราะห์น้ำเสีย. พิมพ์ครั้งที่ 4. กรุงเทพฯ: เรือนแก้วการพิมพ์, 2547.
- กระทรวงอุตสาหกรรม. ประกาศกระทรวงอุตสาหกรรม, พ.ศ. 2566. เรื่อง การจัดการสิ่งปฏิกูลหรือวัสดุที่ไม่ใช้แล้ว. ราชกิจจานุเบกษา. 31 พฤษภาคม 2566. เล่มที่ 140 ตอนพิเศษ 126 ง.
- APHA, AWWA, WEF. Standard Methods for the Examination of Water and Wastewater. 24th ed. Washington, DC: APHA, 2023.
- United States Environmental Protection Agency. Standards of Performance for New Stationary Sources. 40 CFR 60. Appendix A, 2020.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. SW-846, 2014.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Acid Digestion of Sediments, Sludges, and Soils. SW-846 Method 3050B, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Alkaline Digestion for Hexavalent Chromium. SW-846 Method 3060A, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste 3. Physical/Chemical Methods. Separatory Funnel Liquid-Liquid Extraction. SW-846 Method 3510C, 1996.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Ultrasonic Extraction. SW-846 Method 3550C, 2007.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Volatile Organic Compounds in Various Sample Matrices Using Equilibrium Headspace Analysis. SW-846 Method 5021A, 2014.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Purge and Trap for Aqueous Samples. SW-846 Method 5030C, 2003.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Closed System Purge and Trap and Extraction for Volatile Organics in Soil and Waste Sample. SW-846 Method 5035A, 2000.
- United States...

- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Inductively Coupled Plasma-Optical Emission Spectrometry. SW-846 Method 6010D, 2014.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Flame Atomic Absorption Spectrophotometry. SW-846 Method 7000B, 2007.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Arsenic (Atomic Absorption, Gaseous Hydride). SW-846 Method 7061A, 1992.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Chromium, Hexavalent (Colorimetric). SW-846 Method 7196A, 1992.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Mercury in Liquid Waste (Manual Cold Vapor Technique). SW-846 Method 7470A, 1994.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Mercury in Solid or Semisolid Waste (Manual Cold-Vapor Technique). SW-846 Method 7471B, 1998.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Mercury in Solids and Solutions by Thermal Decomposition, Amalgamation, and Atomic Absorption Spectrophotometry. SW-846 Method 7473, 2007.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Selenium (Atomic Absorption, Borohydride Reduction). SW-846 Method 7742, 1994.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Nonhalogenated Organics Using GC/FID. SW-846 Method 8015D, 2003.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Organochlorine Pesticides by Gas Chromatography. SW-846 Method 8081B, 2007.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Polychlorinated Biphenyls (PCBs) by Gas Chromatography. SW-846 Method 8082A, 2007.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Organochlorine Pesticides by Gas Chromatography. SW-846 Method 8081B, 2007.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Chlorinated Herbicides by GC Using Methylation or Pentafluorobenzoylation Derivatization. SW-846 Method 8151A, 1996.
- United States...

- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Volatile Organic Compounds by Gas Chromatography/ Mass Spectrometry. SW-846 Method 8260D, 2018.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Semivolatile Organic Compounds by Gas Chromatography/Mass Spectrometry. SW-846 Method 8270E, 2018.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Cyanide Extraction Procedure for Solids and Oils. SW-846 Method 9013A, 2014.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Cyanide in Waters and Extracts using Titrimetric and Manual Spectrophotometric Procedures. SW-846 Method 9014, 2014.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. pH Electrometric Measurement. SW-846 Method 9040C, 2004.
- United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Soil and Waste pH. SW-846 Method 9045D, 2004.