

ภาคผนวก ง
เอกสารสอบเทียบเครื่องมือ

List of Instrument Certificates for Environmental Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*
1	Analytical Balance	FAT OIL AND GREASE	Mettler Toledo	AB204-S/FACT / 1129361010	United Analyst and Engineering Consultant Co., Ltd.	250422 1 BL002 25	23/4/2025	22/4/2026
			Mettler Toledo	AB204-S/FACT / 1129361010	United Analyst and Engineering Consultant Co., Ltd.	260318-1-BL007-26	3/4/2026	2/4/2027
2	Analytical Balance	TOTAL DISSOLVED SOLIDS	Mettler Toledo	XSR205DU / C210685394	National Food Institute,Ministry of Industry, Thailand	2502226-002-01	20/3/2025	19/3/2026
			Mettler Toledo	XSR205DU / C210685394	United Analyst and Engineering Consultant.Co.,Ltd (UAE)	260311-2-BL001-26	12/3/2026	11/3/2027
3	Analytical Balance	SUSPENDED SOLIDS	Mettler Toledo	XSR205DU / C009071872	National Food Institute,Ministry of Industry, Thailand	2502226-001-01	20/3/2025	19/3/2026
4	BOD Incubator	BIOCHEMICAL OXYGEN DEMAND	ARCO	UC4-1320 / 1021	Technology Promotion Association (Thailand-Japan)	25TM1001	7/7/2025	6/7/2026
5	Continuous Flow Analyzer(CFA)	CYANIDE	Skalar Analytical B.V., the Netherlands	San++5000-02 / 182688	DKSH (Thailand) Ltd.	WO-00074079	23/5/2025	22/5/2026
6	DO Meter	DO	Horiba	LAQUA-DO210 / HE9M0028	Technology Promotion Association (Thailand-Japan)	25TW19	23/1/2025	21/1/2026
7	DO Meter	DO	Horiba	LAQUA-DO210 / HE9M0011	Technology Promotion Association (Thailand-Japan)	26TW102	30/4/2026	30/4/2027
8	DO Meter	BIOCHEMICAL OXYGEN DEMAND	YSI	5100 / 11B 101863	Technology Promotion Association (Thailand-Japan)	25TW29	17/2/2025	16/2/2026
9	DO/BOD Meter	BIOCHEMICAL OXYGEN DEMAND	YSI	4010-2W / 20260326	Technology Promotion Association (Thailand-Japan)v	25TW213	22/10/2025	21/10/2026
10	SCT Meter	SALINITY	Horiba	LAQUA-EC210 / HC9L0013	Technology Promotion Association (Thailand-Japan)	26CH144	2/2/2026	2/2/2027
11	SCT Meter	SALINITY	Horiba	LAQUA-EC210 / HC1L0016	Technology Promotion Association (Thailand-Japan)	25CH168	5/2/2025	3/2/2026

List of Instrument Certificates for Environmental Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*
12	Heating Block	CHEMICAL OXYGEN DEMAND	Hanna Instrument Inc.,Romania	HI839800-02 / 06480019101	Hanna Instruments (Thailand) Ltd.	C26030163	17/3/2026	16/3/2027
13	Heating Block	CHEMICAL OXYGEN DEMAND	Hanna Instruments Italia Srl.	HI 839800-02 / H 018500 I	Hanna Instruments (Thailand) Ltd.	HIT-2510-0375	7/3/2025	6/3/2026
14	Hot Air Oven	SUSPENDED SOLIDS TOTAL DISSOLVED SOLIDS	Memmert	UF55 / B212.0411	Technology Promotion Association (Thailand-Japan)	25TM579	19/3/2025	18/3/2026
15	Inductively Coupled Plasma- Optical Emission Spectrometer(ICP-OES)	LEAD	Agilent Technologies, USA	5110 VDV(G8015AA) / MY18030001	Agilent Technologies (Thailand) Co.,Ltd.	Preventive Maintenance Checklist	3/11/2025	2/11/2026
16	pH Meter and pH Electrode	pH	Mettler Toledo	pH S20 SevenEasyTM / 1231155210	National Food Institute Ministry of Industry, Thailand	2501844-001-01	24/2/2025	23/2/2026
			Mettler Toledo	pH S20 SevenEasyTM / 1231155210	National Food Institute Ministry of Industry, Thailand	2601956-001-01	23/2/2026	22/2/2027
17	UV-VIS Spectrophotometer	NITRATE NITROGEN SULPHIDE	Hitachi	U-2900 / 21E22-009	DQE Services Co.,Ltd.	SP25-021	26/5/2025	25/5/2026
			Hitachi	U-2900 / 21E22-009	DQE Services Co.,Ltd.	SP26-012	18/5/2026	17/5/2027

Due Date of Calibration* : Based on the annual calibration plan. At least 1 time per year.



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3 : EQUIPMENT CALIBRATION AND TESTING SERVICES

534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250

TEL. 0-2717-3000 FAX. 0-2719-9484

Certificate of Testing

Cert.No.: 25TW19

Page.: 1 of 2

Equipment :	DO Meter
Manufacturer :	Horiba
Model :	LAQUA-DO210
Serial No. :	HE9M0028
ID No. :	UAE.EFM.013/2563 (EFM.DO.02/63)
Received Date :	22 January 2025
Test Date :	23 January 2025
Reference :	2501-0774WSC-1
Submitted by :	United Analyst and Engineering Consultant Co.,Ltd. 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260
Laboratory Condition :	Temperature (25 ± 5) °C Humidity (50 ± 20) %
Test Procedure :	In - house method : CP-CH9 by Comparison Technique with Azide Modification Method
Tested by :	Walalak Sirithean 
Approved by :	 Approved Signatory
() Pornthippa Tameyakul () Ponpan Paipim (✓) Saithip Meangmai	
Issue Date :	24 January 2025



Cert.No.: 25TW19

Page.: 2 of 2

Condition of this result of calibration

1. Reference Standard Instruments :

This certification is traceable to the International System of Unit through the reference standards laboratory of Industrial Calibration Center, Technology Promotion Association (Thailand-Japan).

<u>Instruments</u>	<u>Serial No.</u>	<u>ID No.</u>	<u>Certificate No.</u>	<u>Due Date</u>
1. Burette	-	130BU10	23CG1172	22 Mar 2025
2. Balance	14233821	110RC001	24MM131	04 July 2025

2. Standard Material :-

<u>Material</u>	<u>Manufacturer</u>	<u>Lot.No.</u>	<u>Assay</u>
Sodium Thiosulfate 5-Hydrate AR	KEMAUS	2203162447	99.6%

Result : Dissolved Oxygen Meter Adjustment With Air 100 %

Dissolved Oxygen Probe No.: 9K9G0090

Titration Method (Azide Modification Method) (mg/L)	DO Meter Reading (mg/L)	Standard Deviation (mg/L)
8.20	8.21	0.0045

This report was certified only for the instrument we tested. It is allowable to use for study
Intend to use for advertising and referral purpose is prohibited. This report may not be reproduced
other in full, without written approval of the laboratory

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert. No.: 25LM11

Page.: 1 of 2

Equipment : DO Meter with Sensor

Manufacturer : Horiba

Model : LAQUA-DO210

Serial No. : HE9M0028

ID No. : UAE.EFM.013/2563(EFM.DO.02/63)

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong, Bangkok 10260

Location : TPA On Site Calibration Laboratory

Received Order : 22 January 2025

Calibrated Date : 23 January 2025

Ambient Temperature : (26 ± 10) °C

Relative Humidity : (50 ± 30) %

AC Line Voltage : (220 ± 22) V

Calibrated by : Warakorn Lerngagtrakul

Approved by :

Approved Signatory

- () Chakrit Waewwanjua
() Suwit Imjai
(✓) Kunchit Promprat

Issue Date : 29 January 2025

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Equipment : DO Meter with Sensor
Condition As-Received : Used Item
Reference : 2501-0774WSC-2
Procedure Used :-

Cert. No.: 25LM11
Page.: 2 of 2

Calibration were conducted using in-house calibration procedure CP-OT01 according to comparison with Industrial Platinum Resistance Thermometer (IPRT) into Temperature Bath.

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

<u>Instrument</u>	<u>Serial No.</u>	<u>Cert. No.</u>	<u>Traceable</u>	<u>Due Date</u>
1) Digital Thermometer	2188080	2411022	TPA	17 Sep 2025

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This certification is traceable to the International System of Unit.

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function : Temperature measurement.

This instrument was connected with temperature sensor, S/N.: 9K9G0090

<u>Calibration Point</u> (°C)	<u>Immersion Depth</u> (mm)	<u>Standard Temperature</u> (°C)	<u>UUC* Reading</u> (°C)	<u>Error</u> (°C)	<u>Uncertainty</u> (± °C)	<u>Coverage Factor</u> <i>k</i>
15.0	90	15.002	15.0	-0.002	0.16	2.00
30.0	90	30.002	30.0	-0.002	0.16	2.00
45.0	90	45.003	45.0	-0.003	0.16	2.00

UUC* : Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

Calibration Certificate

Certificate No.: 2501844-001-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road,
Bangchack, Prakhonong, Bangkok 10260

Page 1 of 5

Equipment: pH Meter
Manufacturer: METTLER TOLEDO
Model: SevenEasy pH
Serial No.: 1231155210
ID No.: UAE.WAT.010/2553
Order No.: 2501844
Operation No.: 2501844-001
Date of Receipt: 24 February 2025
Date of Calibration: 24 February 2025

Calibrated by Mr.Manas Somsak
Specialist

Approved by 
(Mr.Pheraphat Tuanjit)
Manager, Division of Calibration Laboratory

Date of Issue: 27 February 2025

Responsible for the Technical Management Team

The uncertainties are for a confidence probability of approximately 95%.

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65



Calibration Report

Certificate No.: 2501844-001-01

Equipment: pH Meter **Resolution:** 0.01 pH ; 1 mV

Manufacturer: METTLER TOLEDO **Model:** SevenEasy pH

Serial No.: 1231155210 **Type:** Bench top

ID No.: UAE.WAT.010/2553

Date of Calibration: 24 February 2025 Page 2 of 5

Location: Chemical Calibration Laboratory, National Food Institute

Environment Condition: **Ambient Temperature:** (23.4 ± 1.5) °C **Relative Humidity:** (54 ± 3) %

Condition of Equipment: Good Condition

Condition of this Results of Calibration

1. Calibration Method W-CC-002 : In house method based on direct measurement by using standard voltage calibrator and certified reference material (CRM)

2. Reference Standards / Certified Reference Material

Instruments	Serial / ID No.	Manufacturer	Certificate No.	Due Date
2.1 DC Voltage Calibrator	2709007	Fluke	24E1752	30 May 2025
2.2 Digital Thermometer	2709007	Fluke	2500376-002-01	29 October 2025
2.3 Thermo-Hygro Meter	NFI.BTH 013/23	testo	CC 670420-01	21 May 2025
Certified Reference Material	Lot. No.	Manufacturer	Ref N	Expire Date
2.4 pH buffer 4.008 (Primary pH buffer Solution)	1016435	CPAchem	PH216.L5	25 July 2026
2.5 pH buffer 6.865 (Primary pH buffer Solution)	949186	CPAchem	PH217.L5	30 November 2025
2.6 pH buffer 10.01 (Primary pH buffer Solution)	1016437	CPAchem	PH220.L5	25 July 2025
2.7 pH buffer 7.00 (Standard pH buffer Solution)	C03109	HACH LANGE GmbH	S11M004	16 October 2025

3. This certification is traceable to The International System of Unit (SI Unit)

3.1 Instruments No.2.1	through	NSC-TISI-TIS 17025 Laboratory Accreditation of Calibration No.0008
3.2 Instruments No.2.2 to 2.3	through	NSC-TISI-TIS 17025 Laboratory Accreditation of Calibration No.0061
3.3 Certified Reference Material No.2.4 to 2.6	traceable to	Primary measurement method- Harned cell using calibrated thermometer, barometer, and nanovoltmeter. The Standard Solution preparation and certified by CPAchem Ltd is accredited to ISO 17034 and ISO/IEC 17025
3.4 Certified Reference Material No.2.7	traceable to	PTB Certificate Nr. PTB-PHOA-563/30504/23 and Certificate Nr. PTB-PHOB-555/30620/22 (PTB: Physikalisch-Technische Bundesanstalt, Braunschweig, Germany)

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.



F-CS-012 Revision: 01 Date: 20-04-65



Calibration Report

Certificate No.: 2501844-001-01

Equipment: pH Meter **Resolution:** 0.01 pH ; 1 mV

Manufacturer: METTLER TOLEDO **Model:** SevenEasy pH

Serial No.: 1231155210 **Type:** Bench top

ID No.: UAE.WAT.010/2553

Date of Calibration: 24 February 2025

Page 3 of 5

Calibration Results:

1. Calibration of pH Meter (Manual Temperature Compensation at 25 °C)

Nominal pH	DC Voltage Standard (mV)	Average Indicator Reading		Uncertainty (±mV)	Coverage Factor (k)
		mV	pH		
0	414.122	414	-0.01	0.58	2.00
2	295.815	296	1.99	0.58	2.00
4	177.463	178	4.00	0.58	2.00
6	59.160	59	6.00	0.58	2.00
7	0.001	0	7.00	0.58	2.00
8	-59.159	-59	8.00	0.58	2.00
10	-177.462	-177	10.00	0.58	2.00
12	-295.813	-296	12.00	0.58	2.00
14	-414.121	-414	14.00	0.58	2.00

2. Calibration of pH Meter with Electrode (Manual Temperature Compensation at 25 °C)

Equipment: pH Electrode **Type:** Combined Electrode

Manufacturer: METTLER TOLEDO **Model:** InLab Solids

Serial No.: 3065701 **ID.No.:** N/A

Performance of Electrode system (Three-Point Calibration at pH 4, 7 and 10)

Certified Value @25 °C (pH)	Average Indicator Reading		Relative Slope (%)	Uncertainty (± pH)	Coverage Factor (k)
	pH	mV			
4.008	4.00	165	-	0.0071	2.00
7.001	7.00	-8	97.5	0.0086	2.00
10.010	10.01	-178	95.5	0.0083	2.00
6.876	6.88	0	-	0.0071	2.00



Calibration Report

Certificate No.: 2501844-001-01

Equipment: Digital Thermometer with RTD (pH Meter)

Resolution: 0.1 °C **Model:** SevenEasy pH

Serial No.: 1231155210 **ID No.:** UAE.WAT.010/2553

Manufacturer: METTLER TOLEDO

Date of Calibration: 24 February 2025

Page 4 of 5

Location: Chemical Calibration Laboratory, National Food Institute

Environment Condition:

Ambient Temperature 23.4 °C ± 1.0 °C

Relative Humidity 55.1 % ± 1.7 %

Condition of this results of Calibration:

- Calibration Method :
 - In house method: W-TE-025 by comparison with standard thermometer.
 - The Calibration is determined by comparing with a known temperature from a standard resistance thermometer.
 - The temperature scale in use at this laboratory is the International Temperature scale of 1990 (ITS-90).

2. Reference Standard Instrument :

Instrument	Model	Serial No.	Certificate No.	Due Date	Through
HANDHELD THERMOMETER	1523	2118154	PSL-T 0815/67	24-Jun-25	TISTR
Platinum Resistance Thermometer (PRT)	5627A	877332			

Support Equipment : - Low Temperature Bath (AMETEK RTC-187), Model: RTC-187C , S/N: 670930-00018

- This certificate is traceable to International System of Units (SI Units).
- This certificate was certified only for the instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.

6. Condition of Calibrated item : Good

7. Result of Calibration : ☒ Without adjustment ☐ After adjustment




Calibration Report

Certificate No.: 2501844-001-01

Equipment: Digital Thermometer with RTD (pH Meter)

Resolution: 0.1 °C Model: SevenEasy pH

Serial No.: 1231155210 ID No.: UAE.WAT.010/2553

Manufacturer: METTLER TOLEDO

Date of Calibration: 24 February 2025

Page 5 of 5

Calibration point: 20.0, 25.0 and 30.0 °C

Calibration result:

- The probe was immersed in liquid bath or dry bath to a minimum depth of 120 mm.
- Description of probe, model : N/A S/N : N/A
- Dimension of probe : Diameter 4 mm., Length 120 mm.,
- Sheath material : Stainless Steel

UUC* Reading (°C)	Standard Temperature (°C)	Correction Value (°C)	Uncertainty ± (°C)
20.1	20.001	0.1	0.099
25.1	25.002	0.1	0.099
30.1	30.003	0.1	0.099

Note

- UUC* : Unit Under Calibration

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k= 2, providing a level of confidence of approximately 95 %.

----- End -----

F-CS-012 Revision: 01 Date: 20-04-65



Calibration Certificate

Certificate No.: 2502226-002-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road,
Bangchack, Prakhonong, Bangkok 10260

Page 1 of 4

Equipment: Electronic Balance
Manufacturer: METTLER TOLEDO
Model: XSR205DU
Serial No.: C210685394
ID No.: UAE.WAO.010/2565
Order No.: 2502226
Operation No.: 2502226-002
Date of Receipt: 19 March 2025
Date of Calibration: 20 March 2025

Calibrated by Mr.Yothin Charoensuk
Scientist

Approved by 
(Mr.Pheraphat Tuanjit)
Manager, Division of Calibration Laboratory
Responsible for the Technical Management Team

Date of Issue: 25 March 2025

The uncertainties are for a confidence probability of approximately 95%

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65



Calibration Report

Certificate No.: 2502226-002-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C210685394

ID No.: UAE.WAO.010/2565

Capacity: 82 g / 220 g

Date of Calibration: 20 March 2025

Page 2 of 4

Environment Condition: Ambient Temperature: 21.2 ± 0.6 °C Relative Humidity: 48 ± 3.5 %

Place of Calibration: 208 Balance Room, UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.

Condition of Equipment: Good Condition

Condition of This Results of Calibration:

1. Calibration Method: NFI Method W-MA-001 In-House Method based on UKAS Lab 14 : 2019

2. Reference Standards:

Reference Standard	Model	Serial No.	Calibrated By	Certificate No.	Due Date
Standard Weight Class E2	1mg to 200g	B505567572	TCS	M2404100S	19 April 2025

Instrument	Model	Serial No.	Calibrated By	Certificate No.	Due Date
Thermo-Hygro Meter	608-H1	NFI.BTH 017/23	Quality Reborn	QR25-0542	10 February 2026

3. This certification is traceable to SI UNIT

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

Calibration Results:




1. Repeatability of Reading:

Nominal Value (g)	Standard Deviation of Reading (g)
40	0.0000042
80	0.0000042
100	0.0000000
200	0.0000000

2. Off-Center Error:

A mass of 100 g was placed and moved to various position on pan.

The balance reading obtained is given in the table.

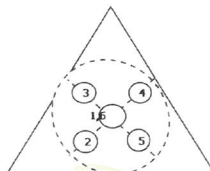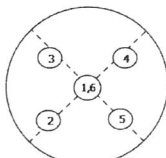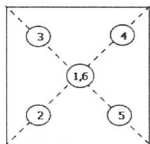

1 (g)	2 (g)	3 (g)	4 (g)	5 (g)	6 (g)	(Maximum Difference) (g)
100.0001	100.0001	100.0001	100.0001	100.0001	100.0001	0.0000

Calibration Report

Certificate No.: 2502226-002-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C210685394

ID No.: UAE.WAO.010/2565

Capacity: 82 g / 220 g

Date of Calibration: 20 March 2025

Page 3 of 4

Calibration Results: (Continued)

Calibration Range: 0-80 g

Calibration Adjustment: Internal Calibration

3. Departure from Nominal Value: (Range: 0 - 82 g ; Resolution: 0.00001 g)

Nominal Value (g)	Standard Value (g)	Average Reading (g)	Correction (g)	Uncertainty (± g)	Coverage Factor k
Unload	0.000000	0.00000	0.00000	0.0000087	2.00
0.001	0.001003	0.00100	0.00000	0.0000090	2.00
0.005	0.005002	0.00501	-0.00001	0.0000092	2.00
0.01	0.010003	0.01002	-0.00002	0.0000089	2.00
0.05	0.049996	0.05001	-0.00001	0.0000096	2.00
0.1	0.100011	0.10002	-0.00001	0.000011	2.00
0.5	0.500016	0.50004	-0.00002	0.000014	2.00
1	1.000003	1.00005	-0.00005	0.000016	2.00
2	2.000023	2.00006	-0.00004	0.000017	2.00
5	5.000015	5.00006	-0.00005	0.000020	2.00
10	10.000009	10.00005	-0.00004	0.000026	2.00
20	20.000030	20.00007	-0.00004	0.000037	2.00
30	30.000039	30.00009	-0.00005	0.000050	2.00
50	50.000028	50.00008	-0.00005	0.000068	2.00
80	80.000067	80.00013	-0.00006	0.00011	2.00



Calibration Report

Certificate No.: 2502226-002-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C210685394

ID No.: UAE.WAO.010/2565

Capacity: 82 g / 220 g

Date of Calibration: 20 March 2025

Page 4 of 4

Calibration Results: (Continued)

Calibration Range: >80-200 g

Calibration Adjustment: Internal Calibration

3. Departure from Nominal Value: (Range: >80 - 200 g ; Resolution: 0.0001 g)

Nominal Value (g)	Standard Value (g)	Average Reading (g)	Correction (g)	Uncertainty (± g)	Coverage Factor <i>k</i>
90	90.00010	90.0002	-0.0001	0.00015	2.00
100	100.00006	100.0001	0.0000	0.00016	2.00
110	110.00007	110.0002	-0.0001	0.00017	2.00
120	120.00009	120.0002	-0.0001	0.00018	2.00
130	130.00010	130.0002	-0.0001	0.00019	2.00
140	140.00013	140.0002	-0.0001	0.00019	2.00
150	150.00009	150.0002	-0.0001	0.00021	2.00
160	160.00010	160.0002	-0.0001	0.00022	2.00
170	170.00012	170.0002	-0.0001	0.00023	2.00
200	200.00013	200.0002	-0.0001	0.00028	2.00

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k* , providing a level of confidence of approximately 95 %.

----- End -----

for N. important

F-CS-012 Revision: 01 Date: 20-04-65



Calibration Certificate

Certificate No.: 2502226-001-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road,
Bangchack, Prakhnong, Bangkok 10260

Page 1 of 4

Equipment: Electronic Balance
Manufacturer: METTLER TOLEDO
Model: XSR205DU
Serial No.: C009071872
ID No.: UAE.WAO.012/2563
Order No.: 2502226
Operation No.: 2502226-001
Date of Receipt: 19 March 2025
Date of Calibration: 20 March 2025

Calibrated by Mr.Yothin Charoensuk
Scientist

Approved by 
(Mr.Pheraphat Tuanjit)
Manager, Division of Calibration Laboratory
Responsible for the Technical Management Team

Date of Issue: 25 March 2025

The uncertainties are for a confidence probability of approximately 95%

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65



Calibration Report

Certificate No.: 2502226-001-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C009071872

ID No.: UAE.WAO.012/2563

Capacity: 82 g / 220 g

Date of Calibration: 20 March 2025

Page 2 of 4

Environment Condition: Ambient Temperature: 21.2 ± 0.6 °C Relative Humidity: 48 ± 3.5 %

Place of Calibration: 208 Balance Room, UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.

Condition of Equipment: Good Condition

Condition of This Results of Calibration:

1. Calibration Method: NFI Method W-MA-001 In-House Method based on UKAS Lab 14 : 2019

2. Reference Standards:

Reference Standard	Model	Serial No.	Calibrated By	Certificate No.	Due Date
Standard Weight Class E2	1mg to 200g	B505567572	TCS	M2404100S	19 April 2025

Instrument	Model	Serial No.	Calibrated By	Certificate No.	Due Date
Thermo-Hygro Meter	608-H1	NFI.BTH 017/23	Quality Reborn	QR25-0542	10 February 2026

3. This certification is traceable to SI UNIT

4. This certificate was certified only for the instrument we calibrated.

5. This result of calibration was found accurate as shown on date and place of calibration only.

Calibration Results:




1. Repeatability of Reading:

Nominal Value (g)	Standard Deviation of Reading (g)
40	0.0000052
80	0.0000042
100	0.0000000
200	0.0000000

2. Off-Center Error:

A mass of 100 g was placed and moved to various position on pan.

The balance reading obtained is given in the table.

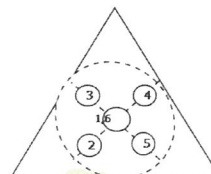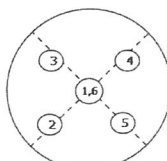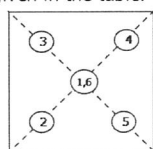

1 (g)	2 (g)	3 (g)	4 (g)	5 (g)	6 (g)	(Maximum Difference) (g)
100.0001	100.0001	100.0001	100.0001	100.0001	100.0002	0.0001

for N. mijkabul

F-CS-012 Revision: 01 Date: 20-04-65



Calibration Report

Certificate No.: 2502226-001-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C009071872

ID No.: UAE.WAO.012/2563

Capacity: 82 g / 220 g

Date of Calibration: 20 March 2025

Page 3 of 4

Calibration Results: (Continued)

Calibration Range: 0-80 g

Calibration Adjustment: Internal Calibration

3. Departure from Nominal Value: (Range: 0 - 82 g ; Resolution: 0.00001 g)

Nominal Value (g)	Standard Value (g)	Average Reading (g)	Correction (g)	Uncertainty (± g)	Coverage Factor k
Unload	0.000000	0.00000	0.00000	0.0000089	2.00
0.001	0.001003	0.00100	0.00000	0.0000092	2.00
0.005	0.005002	0.00500	0.00000	0.0000094	2.00
0.01	0.010003	0.01000	0.00000	0.0000091	2.00
0.05	0.049996	0.05000	0.00000	0.0000098	2.00
0.1	0.100011	0.10000	0.00001	0.000011	2.00
0.5	0.500016	0.50000	0.00002	0.000014	2.00
1	1.000003	1.00001	-0.00001	0.000016	2.00
2	2.000023	2.00005	-0.00003	0.000017	2.00
5	5.000015	5.00005	-0.00003	0.000021	2.00
10	10.000009	10.00005	-0.00004	0.000026	2.00
20	20.000030	20.00012	-0.00009	0.000037	2.00
30	30.000039	30.00012	-0.00008	0.000050	2.00
50	50.000028	50.00014	-0.00011	0.000068	2.00
80	80.000067	80.00020	-0.00013	0.00011	2.00



Calibration Report

Certificate No.: 2502226-001-01

Equipment:

Electronic Balance

Manufacturer: METTLER TOLEDO

Model: XSR205DU

Resolution: 0.00001 g / 0.0001 g

Serial No.: C009071872

ID No.: UAE.WAO.012/2563

Capacity: 82 g / 220 g

Date of Calibration: 20 March 2025

Page 4 of 4

Calibration Results: (Continued)

Calibration Range: >80-200 g

Calibration Adjustment: Internal Calibration

3. Departure from Nominal Value: (Range: >80 - 200 g ; Resolution: 0.0001 g)

Nominal Value (g)	Standard Value (g)	Average Reading (g)	Correction (g)	Uncertainty (± g)	Coverage Factor <i>k</i>
90	90.00010	90.0002	-0.0001	0.00015	2.00
100	100.00006	100.0001	0.0000	0.00016	2.00
110	110.00007	110.0001	0.0000	0.00017	2.00
120	120.00009	120.0002	-0.0001	0.00018	2.00
130	130.00010	130.0002	-0.0001	0.00019	2.00
140	140.00013	140.0002	-0.0001	0.00019	2.00
150	150.00009	150.0002	-0.0001	0.00021	2.00
160	160.00010	160.0002	-0.0001	0.00022	2.00
170	170.00012	170.0002	-0.0001	0.00023	2.00
200	200.00013	200.0002	-0.0001	0.00028	2.00

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k* , providing a level of confidence of approximately 95 %.

----- End -----

F-CS-012 Revision: 01 Date: 20-04-65

for N. ingudat



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3 : EQUIPMENT CALIBRATION AND TESTING SERVICES

534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250

TEL. 0-2717-3000 FAX. 0-2719-9484

Certificate of Testing

Cert.No.: 25TW29

Page.: 1 of 2

Equipment : DO Meter
Manufacturer : YSI
Model : 5100
Serial No. : 11B 101863
ID No. : UAE.WAO.004/2554
Received Date : 14 February 2025
Test Date : 17 February 2025
Reference : 2502-0473DSC-1
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road, Bangchak,
Phrakhanong, Bangkok 10260
Laboratory Condition : Temperature (25 ± 5) °C
Humidity (50 ± 20) %
Test Procedure : In - house method : CP-CH9
by Comparison Technique with Azide Modification Method
Tested by : Walalak Sirithean

Approved by : _____
Approved Signatory
() Chakrit Waewwanjua
() Ponpan Paipim
(✓) Saithip Meangmai
Issue Date : 18 February 2025

เอกสารไม่ควบคุม



Cert.No.: 25TW29

Page.: 2 of 2

Condition of this result of calibration

1. Reference Standard Instruments :

This certification is traceable to the International System of Unit through the reference standards laboratory of Industrial Calibration Center, Technology Promotion Association (Thailand-Japan).

<u>Instruments</u>	<u>Serial No.</u>	<u>ID No.</u>	<u>Certificate No.</u>	<u>Due Date</u>
1. Burette	-	130BU10	23CG1172	22 Mar 2025
2. Balance	14233821	110RC001	24MM131	04 July 2025

2. Standard Material :-

<u>Material</u>	<u>Manufacturer</u>	<u>Lot.No.</u>	<u>Assay</u>
Sodium Thiosulfate 5-Hydrate AR	KEMAUS	2203162447	99.6%

Result : **Dissolved Oxygen Meter Adjustment With Air 100 %**
Dissolved Oxygen Probe No.: 24F100202

Titration Method (Azide Modification Method) (mg/L)	DO Meter Reading (mg/L)	Standard Deviation (mg/L)
8.22	8.22	0.0055

This report was certified only for the instrument we tested. It is allowable to use for study
Intend to use for advertising and referral purpose is prohibited. This report may not be reproduced
other in full, without written approval of the laboratory

-o0o-

เอกสารไม่ควบคุม

เอกสารไม่ควบคุม

บันทึกผลการทวนสอบใบรับรองการสอบเทียบ (Verification of Certificate)							
Certificate No. : 25TW29				Equipment : Do Meter			
Brand : YSI				Model : 5100			
Serial No. : 11B 101863				ID No. : UAE.WAO.004/2554			
Calibration results							
Titration Method	Standart Deviation	Do meter Reading	Error%	Correction%	Error Total Error	Judgement	(Total Error < Judgement)
(mg/L)	(mg/L)	(mg/L)	(mg/L)	(mg/L)	(mg/L)	(± mg/L)	(mg/L)
8.22	0.0055	8.22	0.0000	0.0000	0.0	0.02	pass
ผู้บันทึก.....อิสรา. บุญประกอบ.....ผู้ตรวจสอบ.....ผู้รับรอง.....							
วันที่.....28/02/2025.....วันที่.....28 ก.พ. 68.....							
หมายเลข :							

เก็บในแฟ้ม.....

...../.....

Certificate No.: 250422-1-BL002-25

Code No.: BL002-25

Page: 2 of 3

Equipment: Electronic Balance

Manufacturer: Mettler Toledo

Model: AB204-S/FACT

Readability: 0.0001 g

Serial No.: 1129361010

ID No.: UAE.WAS.002/2552

Max. Capacity: 220 g

Calibration Date: April 23, 2025

Condition As-Received: In Condition

Condition of Equipment:

Condition of This Result of Calibration:

1. Calibration Method: This instrument was calibrated by method UAE.CP.CAL.006 In-House Method based on UKAS Lab 14 : 2022

2. Reference Standards:

Reference Standard:	Model	Serial No.	Calibrated By	Certificate No.	Traceability	Due Date
Standard Weight Class E2 (OIML)	1 mg to 1 kg	B749109122	AMARC	25-009359	Mettler-Toledo	21-Jan-27
Standard Weight Class F1 (OIML)	1 mg to 200 g	11119512	AMARC	24-013840	Mettler-Toledo	04-Feb-26
Instrument	Model	Serial No.	Calibrated By	Certificate No.	Traceability	Due Date
Thermo-Hygro-Baro Meter	MHB-382SD	AK.46457	SUCCESS	SG-H-00997/67	Success Gateway	21-Nov-25
Thermo-Hygro-Baro Meter	MHB-382SD	AK.46457	TPA	25P795	TPA	25-Feb-26

3. This certification is traceable to SI Unit

4. This certification was certified only for the indtrument we calibrated

5. This result of calibration wae found accurate as show on date and place of calibration only.

6. Through the reference standard laboratory of AMARC 25-009359 Calibration 0152

Calibraton Result:




1. Repeatability of Reading:

Nominal Value (g)	Standard Deviation of Reading (g)
200*	0.000045

2. Eccentric or off-center loading

A mass of 100 g was placed and moved to various position on pan

The Balance reading obtained is given in the table.

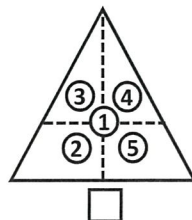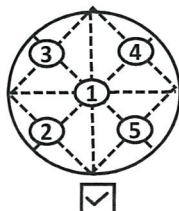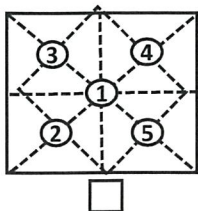

1 (g)	2 (g)	3 (g)	4 (g)	5 (g)	Maximum Difference (g)
100.0000	99.9996	99.9997	100.0003	100.0005	0.0005

เอกสารไม่ควบคุม

Certificate No.: 250422-1-BL002-25

Code No.: BL002-25

Page: 3 of 3

Equipment: Electronic Balance

Model: AB204-S/FACT

Serial No.: 1129361010

Max. Capacity: 220 g

Calibration Date: April 23, 2025

Manufacturer: Mettler Toledo

Readability: 0.0001 g

ID No.: UAE.WAS.002/2552

Calibration Result: (Continued)

Calibration Range: 0 - 200 g

Calibration Adjustment: Internal Calibration

3. Error of indication from nominal or conventional mass value:

Nominal Value (g)	Reference Value (g)	Indication (g)	Correction (g)	Uncertainty ($\pm$ mg)	Coverage Factor k
Unload	0.0000000	0.0000	0.0000	0.10	2.05
0.01	0.0100025	0.0099	0.0001	0.10	2.05
0.05	0.0500056	0.0500	0.0000	0.10	2.05
0.1	0.1000012	0.0999	0.0001	0.10	2.05
0.5	0.5000133	0.5000	0.0000	0.10	2.05
1	1.0000105	1.0000	0.0000	0.10	2.05
10	10.000010	10.0000	0.0000	0.11	2.04
40	40.000076	40.0000	0.0000	0.14	2.00
50	50.000056	50.0000	0.0001	0.13	2.00
80	80.000107	80.0000	0.0001	0.18	2.00
100	100.000109	99.9999	0.0002	0.17	2.00
120	120.00015	119.9999	0.0003	0.21	2.00
150	150.000165	149.9998	0.0003	0.24	2.00
160	160.000175	159.9997	0.0005	0.26	2.00
200	200.000129	199.9998	0.0004	0.30	2.00

4. Effect of Tare test:

Tare Load (g)	Test Load (g)	Indication (g)	Correction (g)
100	20.000041	19.9999	0.0001
	40.000076	39.9998	0.0002
	60.000066	59.9997	0.0003
	80.000107	79.9999	0.0002
	100.000168	100.0004	-0.0003

Remark:

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k, providing a level of confidence of approximately 95%.

เอกสารไม่ควบคุม

o---o-End-o---

Certificate No. : HIT-2510-0375

Page : 1 of 2

CERTIFICATE OF CALIBRATION

Equipment : COD Test Tube Heater

Meter Model : HI839800-02 **Serial No. :** H018500I

Tube Heater : 25 Vial Capacity **Resolution :** 0.1°C

Temperature Range : (-10 to 160)°C **Temperature of Reaction :** 150°C

Manufacturer : Hanna Instruments **Made in :** Romania

Condition As-Received : Used Product **Reference :** RE250401

Ambient Temperature : (25 ± 2)°C **Relative Humidity :** (50 ± 15) % RH

Customer name : United Analyst and Engineering Consultant Co., Ltd.
3 Soi Udomsuk 41, Sukhumvit Rd., Bangchak,
Phrakhanong, Bangkok 10260

Received date : 5 March 2025


Calibrate date : 7 March 2025

Issue date : 7 March 2025

Calibrated Location : Hanna Instruments (Thailand) Ltd.

Calibration Procedure : This calibrator was conducted by using in-house: calibration procedure
CP-04 by using certified reference standard instruments.

Calibrated by : ☒ Mr. Pichit Petthong
☐ Mr. Channarong Soinak

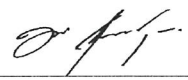Approved by : 
Mr. Anan Suwanchaisakul

Authorized Signatory

 **HANNA**
Instruments
(Thailand) Limited

This certificate was certified only for the instrument we calibrated.

This result of calibration was found accurate on date and place of calibration only.

** This certificate may not be reproduced other than in full, except with the prior written **

approval of the head of Hanna Instrument (Thailand) **เอกสารไม่ควบคุม**

Condition of this calibration result:

Reference Standard Instruments : This certification is traceable to the international unit of unit maintained through:

Instruments	Model	Serial No.	Certificate No.	Traceable
Data Acquisition Switch Unit	34970A	MY44065265	WK2407-141-1	WK Electric Co., Ltd.
Digital Thermo-Hygrometer	HT-771SD	AI.07155	25H171	Technology Promotion Association (Thailand-Japan).

Calibration Result:

Measurement Temperature Source Accuracy for COD Reactor.

Capacity (Vial)	Nominal Value (°C)	Average Value (°C)	Uncertainty of Measurement (±°C)
25 Vial	150.0	150.4	0.47

Unit : °C

(1A)	(2A)	(3A)	(4A)	(5A)
150.407	150.377	150.269	150.402	150.422
(1B)	(2B)	(3B)	(4B)	(5B)
150.426	150.394	150.644	150.690	150.542
(1C)	(2C)	(3C)	(4C)	(5C)
150.477	150.303	150.627	150.257	150.176
(1D)	(2D)	(3D)	(4D)	(5D)
150.462	150.456	150.199	150.406	150.102
(1E)	(2E)	(3E)	(4E)	(5E)
150.185	150.513	150.235	150.460	150.442

Figure: Shows the location of the temperature source.

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k = 2$, providing a level of confidence of approximately 95%

** End of certificate **

เอกสารไม่ควบคุม

Service Report

TO	FOR
Company: United Analyst and Engineering Consultant Co., Ltd. _ Bangkok-HQ Address: 700/2 หมู่ที่ 1 Phrakhanong District, Bangkok, 10260	Work Order Number: WO-00074079 Contact: Karnphong Boonpuang Email: karnphong.b@uaeconsultant.co.th Tel: +66 2763 2828 (7021), +66 8 6347 7390

WORK ORDER INFORMATION			
Top-Level		Order Type	Preventive Maintenance
Installed Product ID	IB-00105024	Billing Type	Chargeable
Product	SKALAR 2SAN59000 SAN++ Classic CFA 230V 2SAN59000	PO No.	HPO-250400209
Serial No.	182688	Warranty No.	
		Contract No.	

PRODUCTS SERVICED		
Installed Product Id	Serial Number	Product
IB-00105024	182688	SKALAR 2SAN59000 SAN++ Classic CFA 230V 2SAN59000

PROBLEM DESCRIPTION
PM 1 ครั้ง/ปี **ใบเสนอราคาเลขที่ Q-120095

Line Number	Engineer	Start Date And Time	End Date And Time	Billable Labor Hour	Billable Travel Hour	Travel KM
WL-00342192	Yongyuth Chanphong	05/23/2025 9:30 AM	05/23/2025 6:00 PM	8.5		
Total				8.5	0	0

Line Number	Work Description
WL-00342192	- ทำ PM เสร็จแล้ว, เครื่องพร้อมใช้งาน

[illegible]

EXPENSES			
Part No	Expense Type	Description	Line Quantity

RECOMMENDED PARTS
อะไหล่สำรองที่ควรสั่งซื้อทั้งหมด 7 รายการ คือ Pump tube 3 รายการ (SA3028, SA3032 และ SA3034), หลอดไฟ Halogen 6V/10W (90020012) 1 รายการ จำนวน 2 หลอด, Tubing polyethylene 3 รายการ (SA3142, SA5141 และ SA5142)

REMARKS

Travel Time Disclaimer:

Please note that the travel time in this report only includes time taken to reach the installed equipment location. It does not include our engineer's return travel time.

Customer Signature:

K. Ranphong

Technician: Yongyuth Chanphong

Job Title: Service Manager

Email: yongyuth.yc@dksh.com

Customer Signature

Date: 06/06/2025

Test Report

Customers	United Analyst and Engineering Consultant Co., Ltd.		
Equipment	Continuous Flow Analyzer	Manufacturer	SKALAR
Controller Mdel	SA5000	Auto Sample Model	SA1052
Controller Serial No.	182688	Auto Sample Serial No.	181729
Date of test	23-May-2025	Period	12 Month
Environment temperature	24.3 °C	Humidity	54.4 %RH

Results

Instrument Checked

Instrument Checked

Item	Characteristic	Before		After		Remark	
1	Visual inspect	☒ Pass	☐ Fail	☒ Pass	☐ Fail		
2	Power supply (210 - 240 VAC)	220	VAC	220	VAC		
3	Computer	☒ Pass	☐ Fail	☒ Pass	☐ Fail		
4	Program	☒ Pass	☐ Fail	☒ Pass	☐ Fail		
5	Auto sampler	☒ Pass	☐ Fail	☒ Pass	☐ Fail		
6	Module holder						
	- Motor pump	☒ Pass	☐ Fail	☒ Pass	☐ Fail		
	- Pump tube	☐ Pass	☒ Fail	☒ Pass	☐ Fail	*	
	- Air-injection	☒ Pass	☐ Fail	☒ Pass	☐ Fail	*	
	- Chemistry manifolds, Switching valve, Coil, Membrane	☒ Pass	☐ Fail	☒ Pass	☐ Fail		
7	Detector						
	- Filter	☒ Pass	☐ Fail	☒ Pass	☐ Fail		
	- Flow cell	☒ Pass	☐ Fail	☒ Pass	☐ Fail		
	- Lamp	☒ Pass	☐ Fail	☒ Pass	☐ Fail		
8	Interface	☒ Pass	☐ Fail	☒ Pass	☐ Fail		
9	Rinsing valves	☐ Pass	☐ Fail	☒ N/A	☐ Pass	☐ Fail	☒ N/A
10	Temperature / Reactor	☒ Pass	☐ Fail	☐ N/A	☒ Pass	☐ Fail	☐ N/A
11	Flame photometer	☐ Pass	☐ Fail	☒ N/A	☐ Pass	☐ Fail	☒ N/A
12	UPS / Stabilizer	☒ Pass	☐ Fail	☐ N/A	☒ Pass	☐ Fail	☐ N/A

Warning and Error Checked

Item	Event	Before		After	
13	Error list	☒ None		☒ None	
		☐ Appear : _____		☐ Appear : _____	

Check with Standard

Item	Characteristic	Before	After	Remark
14	Base Line Test	☐ Pass ☐ Fail ☒ N/A	☒ Pass ☐ Fail ☐ N/A	
15	Detector Signal Test	☐ Pass ☐ Fail ☒ N/A	☒ Pass ☐ Fail ☐ N/A	

Summary of checked

- ☒ The instrument can work normally and efficiently. (เครื่องมือวัดสามารถทำงานได้ปกติและมีประสิทธิภาพ)
- ☐ The instrument can work but it's requiring to maintenance. (เครื่องมือวัดสามารถทำงานได้แต่ต้องบำรุงรักษา)
- ☐ The instrument could not work it's requiring to repair. (เครื่องมือวัดไม่สามารถทำงานได้แต่ต้องการซ่อมบำรุง)

Remark :

* Pump tube, Tubing polyethylene และ Air tube เริ่มเสื่อมสภาพ ได้เปลี่ยนอะไหล่ทั้งหมดแล้วตามระยะการใช้งาน

หมายเหตุ แนะนำอะไหล่ที่ควรสั่งเพิ่มเติม 7 รายการ ดังนี้

1. อะไหล่ พารามิเตอร์ Ammonia จำนวน 3 รายการ (SA3032, SA5141 และ 90020012)
2. อะไหล่ พารามิเตอร์ Phenol และ Cyanide จำนวน 6 รายการ (SA3028, SA3034, SA3142, SA5142 และ 90020012)

Standard Equipment Used

Equipment	Equipment I.D.	
Digital multi meter	S/N 57600592	Due date : 19-Jun-2025
Thermo hygrometer	S/N 39520444/904	Due date : 27-Dec-2025

Test By :

(Mr. Yongyuth Chanphong)

Approved by :


(Mr. Eknamong Wankliang)



Position :



Supervisor, Technical Service



Position :

Manager, Technical Services

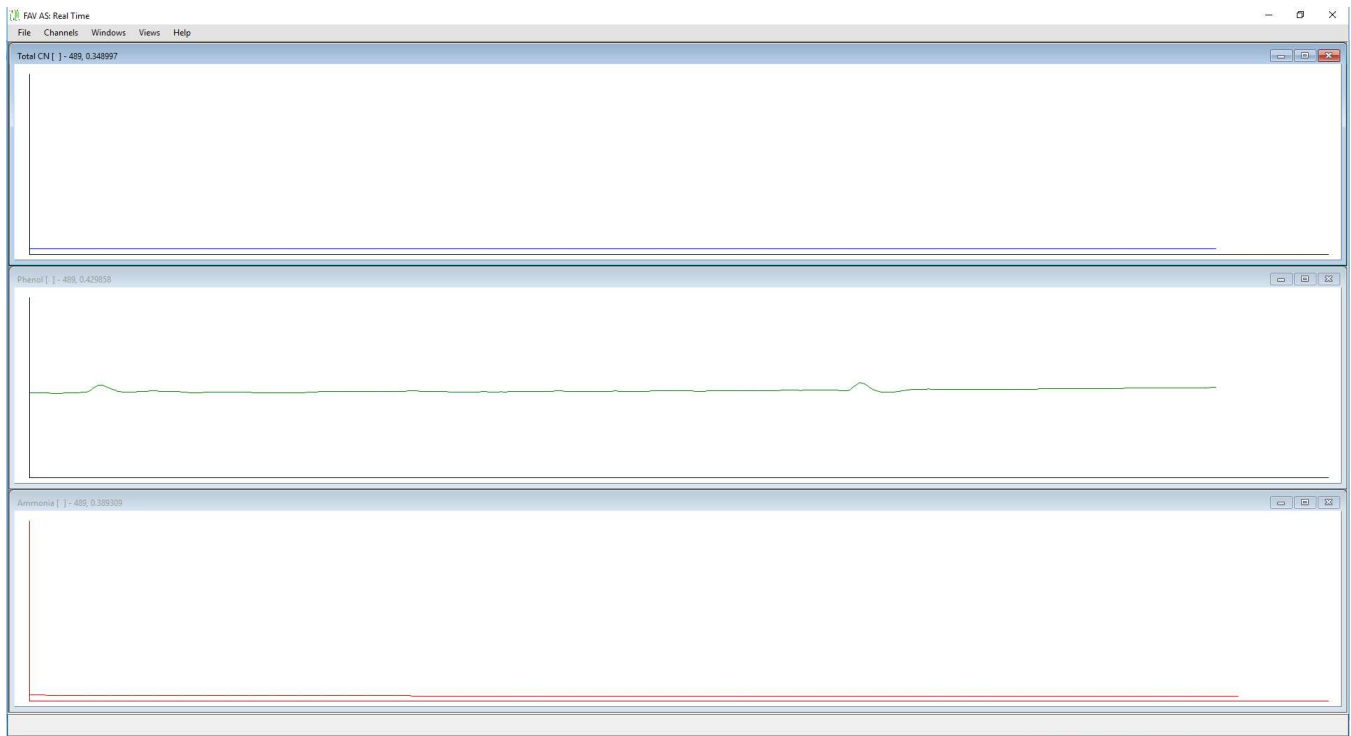Base Line Test : Reagent_Baseline_CN_Phenol_NH3

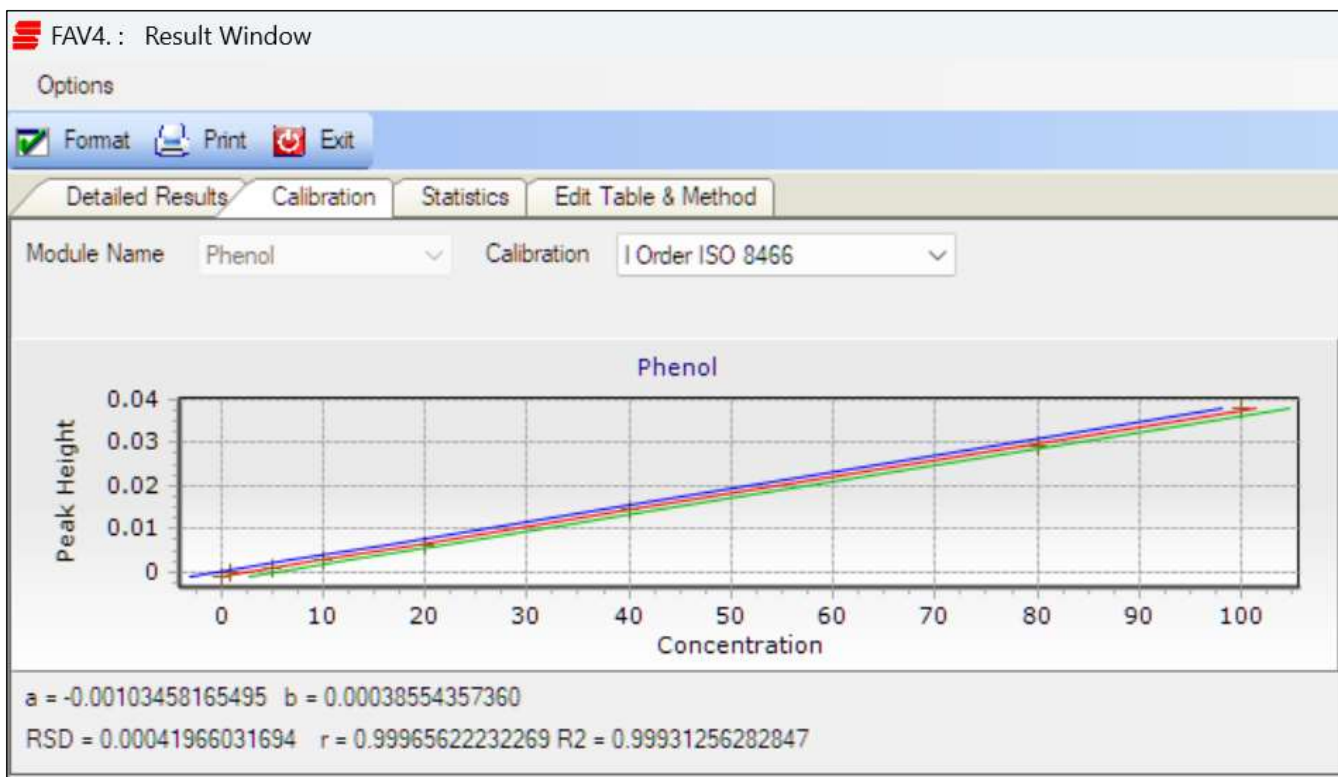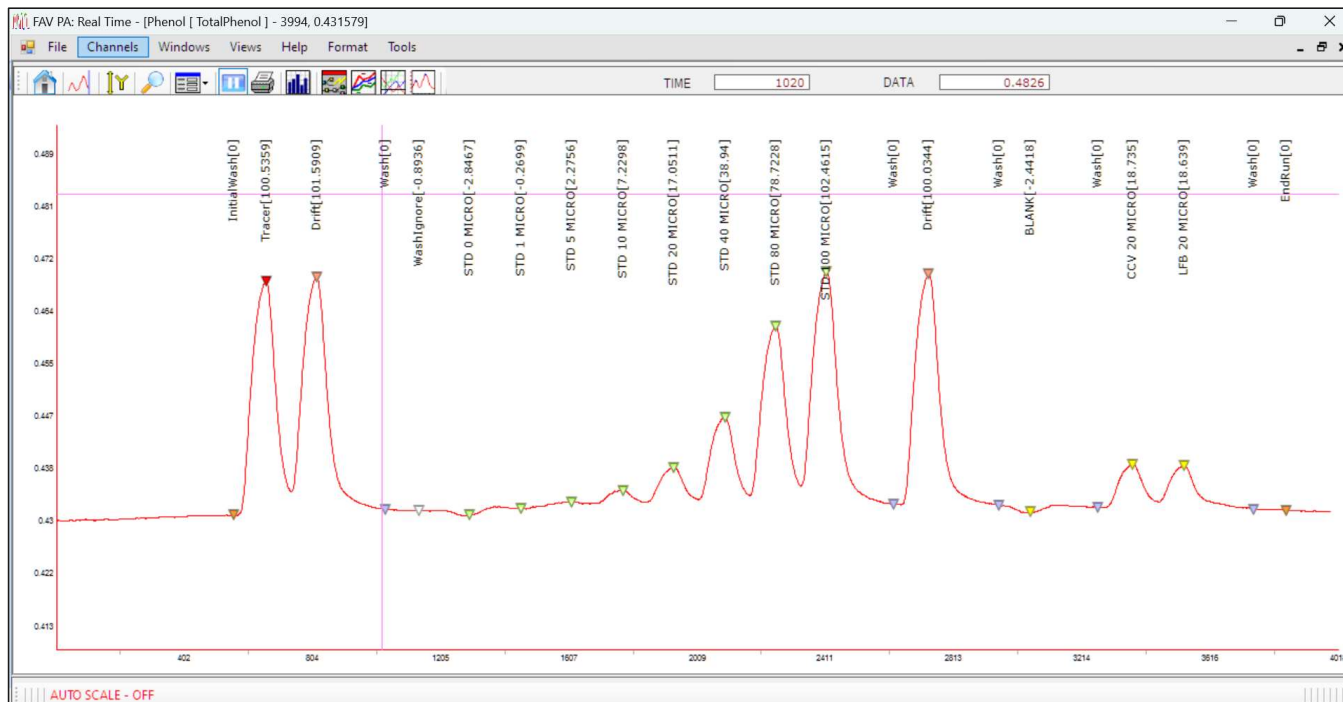Detector Signal Test : Phenol

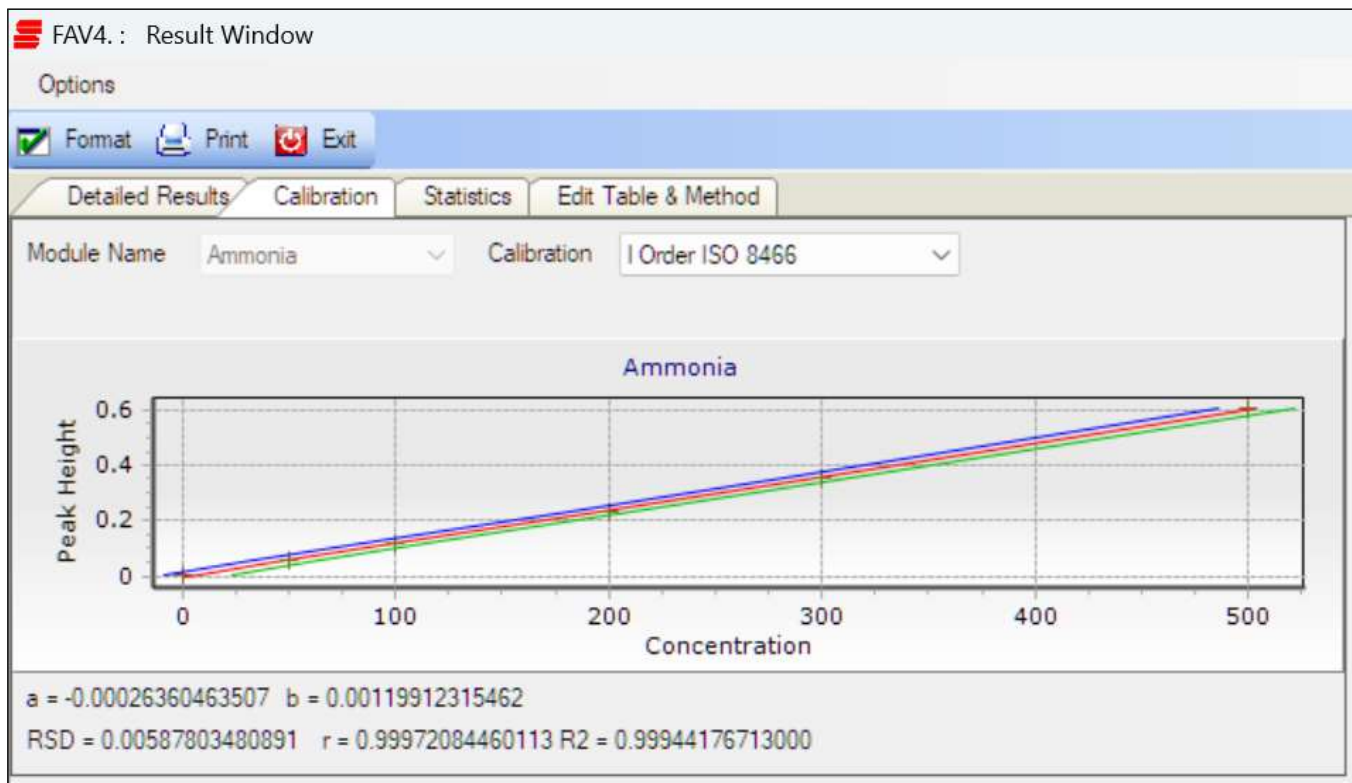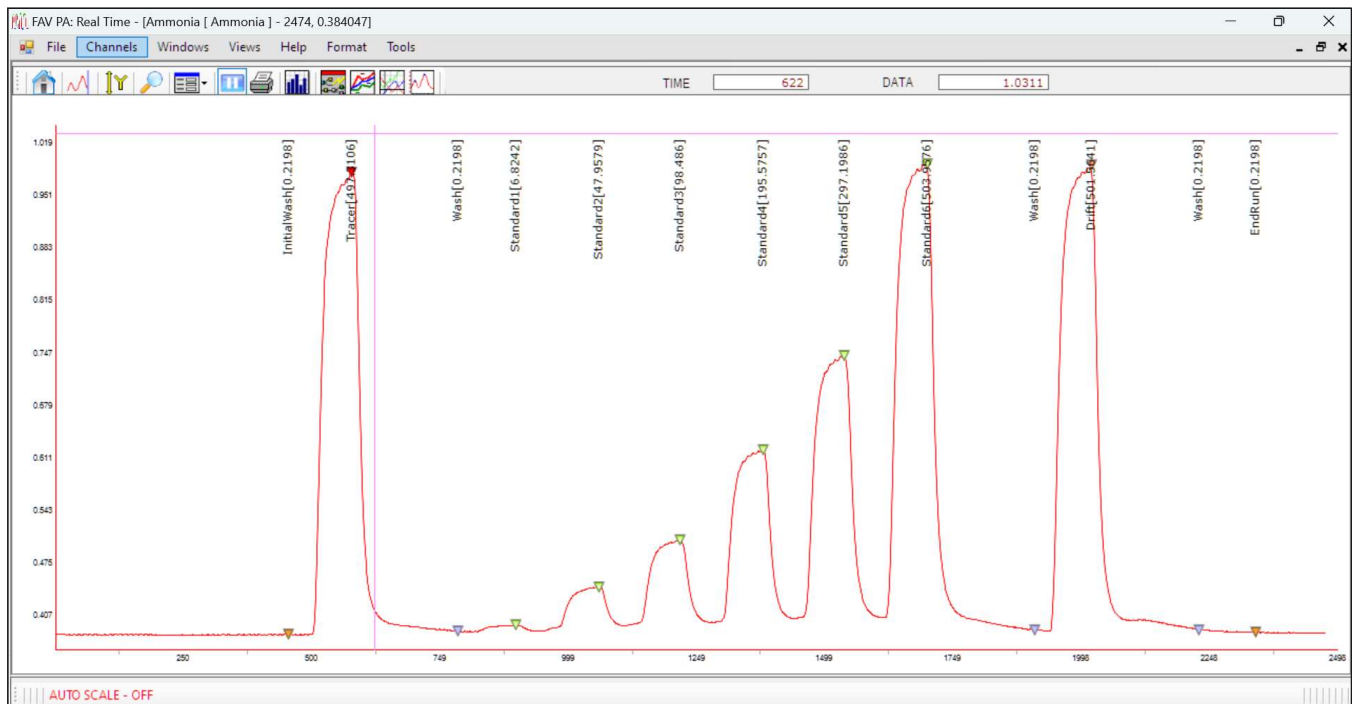
Detector Signal Test : NH3

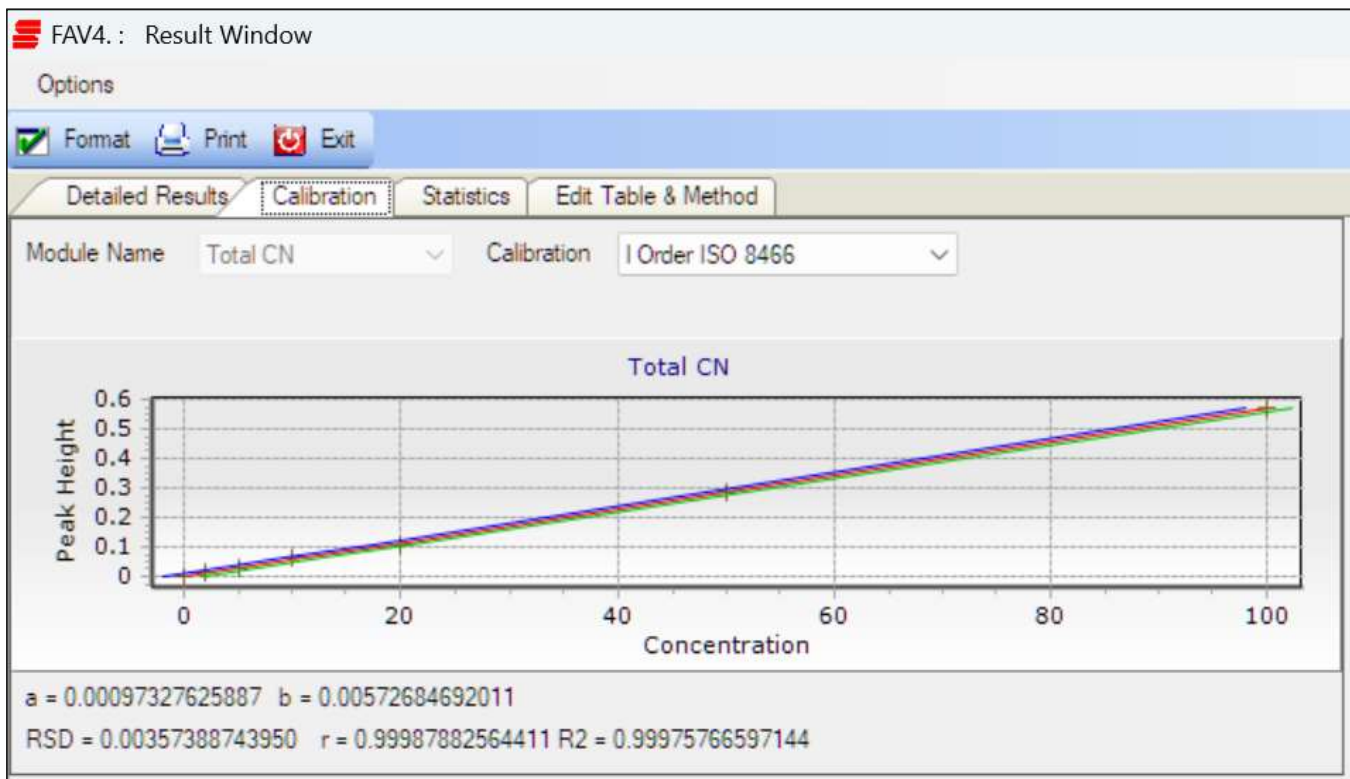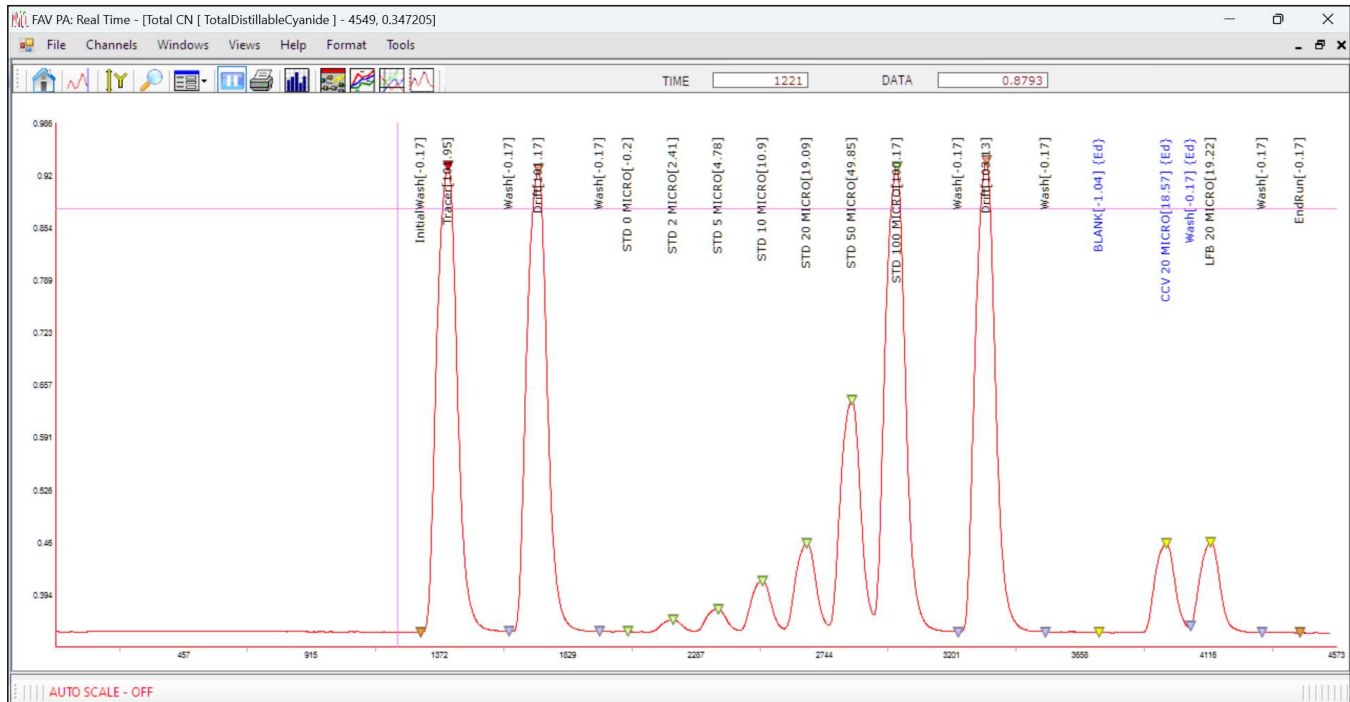
Detector Signal Test : CN




TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert. No.: 25TM579

Page : 1 of 3

Equipment : Hot Air Oven

Manufacturer : Memmert

Model : UF 55

Serial No. : B212.0411

ID No. : UAE.WAO.005/2556

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260

Location : Lab Floor 2

Received Order : 19 March 2025

Calibration Date : 19 March 2025

Ambient Temperature : (26 ± 10) °C

Relative Humidity : (50 ± 30) %

AC Line Voltage : (220 ± 22) V

Calibrated by : Man Pattanapongpaiboon

Kunchit

Approved by :

Approved Signatory

- () Chakrit Waewwanjua
() Suwit Imjai
(✓) Kunchit Promprat

Issue Date : 27 March 2025

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : Hot Air Oven
Condition As-Received : Used Item
Reference : 2503-04370C-3

Cert. No.: 25TM579

Page : 2 of 3

Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD) and Thermocouple Type T.

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

<u>Instrument</u>	<u>Serial No.</u>	<u>Cert. No.</u>	<u>Traceable</u>	<u>Due Date</u>
1) Data Acquisition	MY44073381	24LM73	TPA	18 May 2025

2. This certificate is valid only to the item calibrated on date and place of calibration.


3. This certification is traceable to the International System of Unit.

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Close

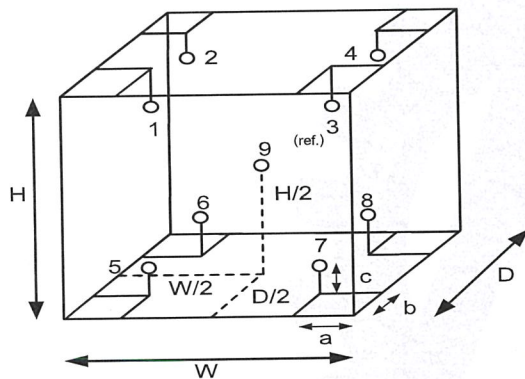

Environment during calibration		
	Beginning	Finished
Temp. (°C)	27	28
REL.Humid. (%)	49	55
AC Supply (Volt)	221	224

Ref. Std. ID No.: @ Calibration Point		
Position :	(120,180) °C	(104) °C
1	23-01TC-01	1RTD-2/1
2	23-01TC-02	1RTD-2/2
3	23-01TC-03	22-01RTD-03
4	23-01TC-04	1RTD-2/4
5	23-01TC-05	1RTD-2/5
6	23-01TC-06	1RTD-2/6
7	23-01TC-07	23-01RTD-07
8	23-01TC-08	1RTD-2/8
9 (ref.)	23-01TC-09	23-01RTD-09

Probe Installation Details :

a = 5.0 cm
b = 5.0 cm
c = 5.0 cm

Dimension of Chamber :

D = 0.50 m
W = 0.80 m
H = 0.75 m
Capacity = 0.30 m³



Equipment : Hot Air Oven
Condition As-Received : Used Item
Reference : 2503-0437OC-3
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Close

Cert. No.: 25TM579

Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor <i>k</i>
104.0	104.0	104.0	0.040	0.43	0.78	2
120.0	120.0	120.0	0.64	1.3	1.6	2
180.0	180.0	180.0	0.49	1.5	1.8	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (±°C)
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	
104.0	104.335	104.135	104.363	104.317	103.649	103.738	104.179	104.229	104.025	0.42
120.0	119.575	119.366	119.807	119.905	118.994	119.194	119.888	119.994	120.064	1.1
180.0	180.286	179.510	180.401	180.551	179.281	179.463	180.196	180.451	180.374	1.2

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-



Certificate of Calibration

Cert.No.: 25CH168

Page.: 1 of 3

Equipment : Conductivity Meter
Manufacturer : Horiba
Model : LAQUA-EC210
Serial No. : HC1L0016
ID No. : UAE.EFM.014/2565(EFM.SCT.01/65)
Condition As-Received: Used Item
Received Date : 04 February 2025
Calibration Date : 05 February 2025
Reference : 2502-0107WSC-4
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong, Bangkok 10260
Ambient Temperature : $(25 \pm 2.5) ^\circ\text{C}$
Relative Humidity : $(50 \pm 15) \%$
Calibration Procedure: In -house method :
- CP-CH6 by direct measurement
with certified reference material (CRM)
- CP-CH8 by comparison with temperature standard
Calibrated by : Warakorn Lerngagtrakul
Approved by : 
Approved Signatory
() Chakrit Waewwanjua
() Ponpan Paipim
(✓) Saithip Meangmai
Issue Date : 06 February 2025

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Cert.No.: 25CH168

Page.: 2 of 3

Condition of this result of calibration

1. Reference Standard Instrument :-

<u>Instrument</u>	<u>Serial No.</u>	<u>ID No.</u>	<u>Certificate No.</u>	<u>Due date</u>
1) Thermometer	1963878	130RC095	24I995	09 Sep 2025
2) Ref. Std.Thermometer	4982054	110RC044	24I757	14 July 2025

- This Certification is traceable to SI Through Technology Promotion Association (Thailand - Japan)

2. Certified Reference Materials :-

- Conductivity calibration solution, CPA chem Ltd., The measurement results are traceable to SI through CPA chem Ltd., ANSI-ASQ National Accreditation Board, Accredited No. AR-1835

<u>Conductivity Solution</u>	<u>Manufacturer</u>	<u>Lot No.</u>	<u>Exp. date</u>
1412.9 μ S/cm	CPA Chem	1005307	15 June 2025
12.881 mS/cm	CPA Chem	1005308	15 June 2025

- Control Conductivity calibration solution temperature by Water bath (25 ± 0.1) °C

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration results

Function : Conductivity Measurement

(*) After Adjustment at 1412.9 μ S/cm

Conductivity Electrode Serial No.: 9B1J0042

Standard Conductivity Solution	Before Adjustment UUC* Reading	After Adjustment UUC* Reading	Uncertainty of Measurement ($\pm$)	Coverage factor <i>k</i>
1412.9 μ S/cm	1355 μ S/cm	1413 μ S/cm	9.2 μ S/cm	2.00
12.881 mS/cm	12.26 mS/cm	12.72 mS/cm	0.086 mS/cm	2.00

Remark : - UUC* = Unit Under Calibration

เอกสารไม่ควบคุม



Cert.No.: 25CH168

Page.: 3 of 3

Calibration Results

Function : Temperature Measurement

This equipment was connected with Temperature Probe;

- Model : 9383
- Serial No. : 9B1J0042

Dimension of probe;

- Length : 110 mm
- Diameter : 16 mm
- Immersion Depth : 100 mm

Calibration Result : Without adjustment

Calibration Point (°C)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty of Measurement (± °C)	Coverage factor <i>k</i>
15.0	15.003	15.0	-0.003	0.13	2.00
30.0	30.002	29.9	-0.102	0.13	2.00
45.0	45.001	44.9	-0.102	0.13	2.00

Remark : - UUC* = Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert. No.: 25TM1001

Page : 1 of 3

Equipment : BOD Incubator
Manufacturer : ARCO
Model : UC4-1320
Serial No. : -
ID No. : UAE.WAO.018/2559
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260
Location : Lab Floor 2
Received Order : 07 July 2025
Calibration Date : 07 July 2025
Ambient Temperature : (26 ± 10) °C
Relative Humidity : (50 ± 30) %
AC Line Voltage : (220 ± 22) V

Calibrated by : Man Pattanapongpaiboon

Kunchit

Approved by :

Approved Signatory

- () Chakrit Waewwanjua
() Suwit Imjai
(✓) Kunchit Promprat

Issue Date : 17 July 2025

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2507-0146OC-2

Cert. No.: 25TM1001

Page : 2 of 3

Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD).

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

<u>Instrument</u>	<u>Serial No.</u>	<u>Cert. No.</u>	<u>Traceable</u>	<u>Due Date</u>
1) Data Acquisition	MY59003411	24LM192	TPA	24 Dec 2025

2. This certificate is valid only to the item calibrated on date and place of calibration.


3. This measurement result is traceable to the International System of Unit maintained through :

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Not Available

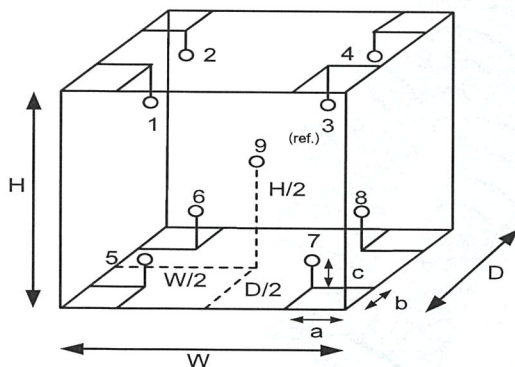

Environment during calibration		
	Beginning	Finished
Temp. (°C)	28	28
REL.Humid. (%)	54	57
AC Supply (Volt)	224	226

Position :	Ref. Std. ID No.:
1	25-20RTD-2/1
2	25-20RTD-2/2
3	25-20RTD-2/3
4	20RTD-2/4
5	20RTD-2/5
6	20RTD-2/6
7	20RTD-2/7
8	20RTD-2/8
9 (ref.)	20RTD-2/9

Probe Installation Details :

a = 10 cm
b = 10 cm
c = 10 cm

Dimension of Chamber :

D = 0.62 m
W = 1.2 m
H = 1.2 m
Capacity = 0.89 m³



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2507-0146OC-2
Result of Calibration :- (*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Not Available

Cert. No.: 25TM1001

Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor <i>k</i>
20.0	20.0	19.9	0.18	0.60	0.83	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	
20.0	20.202	19.900	20.390	20.077	20.111	20.044	20.275	20.173	19.835	0.38

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

Agilent CrossLab Start Up Services

Agilent 5100 5110 ICP-OES Preventive Maintenance

Agilent Preventive Maintenance provides factory recommended service for your analytical instruments to assure reliable operation and the accuracy of your results

Delivered by highly trained and certified service engineers using genuine Agilent parts and supplies, Agilent Preventive Maintenance provides what you need to reduce unplanned downtime and keep your systems operating at their peak performance.

This checklist is used as a guide for completing the preventive maintenance tasks. A signed copy of this checklist is provided for your records.

Introduction

Customer Information

- Customers should provide all necessary operating supplies upon request of the engineer.
- A customer representative should be available to the engineer while performing the preventive maintenance procedures. Customers are responsible for regular maintenance and are encouraged to observe the service representative.
- Any parts not included in the Parts Lists section of this document are not part of the recommended Preventive Maintenance service nor are they included in the price of this service.
- If a system requires the use of extra or special procedures and/or parts for the maintenance service, then these must be ordered separately and charged as a repair, which may incur additional costs.
- For customers using HF applications, the instrument should be returned to its standard sample introduction system.

Important Customer Web Links

- To access **Agilent University**, visit <http://www.agilent.com/crosslab/university/> to learn about training options, which include online, classroom and onsite delivery. A training specialist can work directly with you to help determine your best options.
- To access the **Agilent Resource Center** web page, visit <https://www.agilent.com/en-us/agilentresources>. The following information topics are available:
 - ▶ Sample Prep and Containment
 - ▶ Chemical Standards
 - ▶ Analysis
 - ▶ Service and Support
 - ▶ Application Workflows
- The **Agilent Community** is an excellent place to get answers, collaborate with others about applications and Agilent products, and find in-depth documents and videos relevant to Agilent technologies. Visit <https://community.agilent.com/welcome>
- Videos about specific preparation requirements for your instrument can be found by searching the **Agilent YouTube** channel at <https://www.youtube.com/user/agilent>
- **Need to place a service call?** Flexible Repair Options | Agilent

Service Engineer's Responsibilities

- Contact the customer and ensure that all necessary supplies are available before the preventive maintenance visit.
- Only select those pages that relate to the system or module being serviced.
- Complete empty fields with the relevant information.
- Complete the relevant checkboxes in the checklist using either a "X" or tick mark "✓".
- Check **"Service not applicable"** check boxes to indicate services/tasks not delivered, as appropriate.
- Complete the Preventive Maintenance services in the most logical order relevant to the individual system service in the order of the tasks listed.
- Complete the **Service Review** section together with the customer.
- Complete the fields for page numbers at the foot of each selected page
- Add relevant page numbers to selected pages and complete the total number of pages field in the Service Completion section
- **Ask the customer to sign the Service Verification section including the customer's and your signature.**

Instrument Maintenance

System Information

☐ Check this box if an instrument configuration report is attached instead of completing the table.

Instrument System Name and ID
Instrument System Site and Location

5110 VDV ICP-OES
UAE Consultant

List System Component Product Numbers	List the Serial Numbers of each Component
1. G 8015 A	MY 18030001
2.	
3.	
4.	
5.	
6.	
7.	
8.	
9.	

ICP-OES Configuration Table	Circle the type or write in the type if other
Nebulizer Type	SeaSpray OneNeb Conikal Other
Spray Chamber	Cyclonic Single Pass Cyclonic Double Pass Other
Torch	Radial Dual View Other
Torch Type	One Piece Semi Demountable Fully Demountable Other
Injector Diameter	2.4mm 1.8mm 1.4mm 0.8mm Other
Injector Material	Quartz Ceramic Other

Preparation

- ☒ Discuss any specific issues with the customer before starting.
- ☒ Review the instrument logbook for recorded problems and comments.
- ☒ Save instrument control settings before starting the procedure.
- ☒ Perform a general inspection of the system for cleanliness.
- ☒ Check for proper installation of parts, assemblies, sensors etc.
- ☒ Check system for required installation of components and implementation of Service Notes
- ☒ Check for required firmware/software updates and verify with customers if they would like them installed.
- ☐ For HF application systems, if standard sample introduction system was not installed, ask the customer to install it. *N/A*
- ☒ Ask the customer to remove any samples from the ICP-OES sample introduction area, auto sampler or around the ICP-OES.

Preventive Maintenance Procedures

Record Pre-PM instrument performance

- ☒ Run Instrument Performance test.
- ☒ Record results in Instrument Performance Test Results Table – Pre-PM.

Clean and inspect ICP-OES system

- ☒ Look for any obvious external damage or problems.
- ☒ Inspect water cooling hoses, gas lines and power cord for excessive wear or damage.
- ☒ Perform a general internal inspection of the system for excessive dust accumulation, clean if necessary.
- ☒ Inspect sample introduction components and record any required maintenance in the Service Engineer Comments and notify the customer as the required actions required.
- ☒ Record the instrument operating conditions in the ICP-OES Status Results Table.
- ☒ Replace the polychromator purge filter.
- ☒ Replace the radial pre-optics window
- ☒ Replace the axial pre-optics window for SVDV and VDV instruments.
- ☒ Check exhaust flow for the correct positive extraction at the exhaust duct to insure they meet minimum specifications.
- ☒ Replace air inlet dust filter.
- ☐ Replace high capacity air inlet dust filter element if installed. N/A
- ☒ Remove and clean instrument water inlet filter.

Agilent Water Recirculator

- ☐ **Service not applicable**
- ☒ Drain cooling fluid and remove any particles from the chiller reservoir
- ☒ Remove, clean and reinstall water inlet metal mesh filter if present.
- ☒ Re fill with Agilent Cool Clear cooling fluid.
- ☒ Clean the cooling system Air filter and the condenser.

SPS 3 Auto Sampler

- ☒ **Service not applicable**
- ☐ Power cycle the autosampler and verify successful initialization.
- ☐ Inspect X and Z axis belts for wear. Replace is necessary.
- ☐ Clean X and Z axis slide shafts.
- ☐ Using customer's racks and the Agilent software move the sample probe to the 4 outermost corners and rinse port, ensure that the probe is approximately centered in the vial.

SPS 4 Auto sampler

- ☒ **Service not applicable**
- ☐ Clean the spill tray, rack location mat, end frames and chassis with a damp soft cloth and diluted mild detergent.
- ☐ Clean the auto sampler cover panels, if cover kit is installed, with domestic window cleaner.
- ☐ Check the X-axis and Z-axis drive belts for cracks, splits, damaged teeth, excessive fraying, color changes or degradation from fumes.
- ☐ Check the X-axis, Theta-axis and Z-axis FFC cables for cracks, incorrect positioning, damaged edges or damaged connectors.
- ☐ Pump Tubing Replacement. Replace peristaltic pump tubing. Replace all tubing that goes from the rinse station to the pump and from the pump to the waste/rinse bottles
- ☐ Test using customer's tray and move the sample probe to the sample vial 1, wash vial and rinse port and ensure that the probe is centered in the vial. If not use calibration wizard and calibrate the position.

AVS 4, 6, 7 Advanced Valve System

- ☒ **Service not applicable**
- ☐ Replace valve rotor seal
- ☐ Check fittings for signs of leaks
- ☐ Check tubing including autosampler tubing for kinks or excessive wear
- ☐ Check high flow pump for signs of leaks

ICP-OES adjustment

- ☒ Check position of Zn peak, adjust if required.
- ☒ Check Argon Ratio, adjust to specified value if required.
- ☒ Perform Detector Calibration.
- ☒ Perform Instrument Calibration.

Record Post-PM instrument performance

- ☒ Run Instrument Performance test.
- ☒ Record results in Instrument Performance Test Results Table - Post PM.
- ☒ For systems using ICP Expert version 7.3 and above, run the following Instrument tests
 - ☒ Subsystem Communications Test
 - ☒ Air Flow
 - ☒ Water Flow
 - ☒ Gas Flows
 - ☒ RF Generator
 - ☒ Camera Test
 - ☒ Optics Test
 - ☒ Nebulizer Test
- ☒ Record the result in the Instrument Test Results Table

Restore Instrument

- ☐ For HF applications, ask the customer to reinstall their sample introduction system. *N/A*
- ☒ Leave system in an idle state: on and purging.
- ☒ Guidance: If the PM service is performed prior to a qualification service, then use the qualification procedure as a guide for final instrument set up and checkout.

Service Review

- ☒ Attach available reports/printouts of all tests to this documentation.
- ☒ Record the Preventive Maintenance service activity in the customer's records/logbook.
- ☒ Record the PM event in the Smart Alerts logbook, if applicable.
- ☒ Update/reset instrument maintenance counters as appropriate.
- ☒ Affix the PM sticker to the system or instrument logbook based on the customer's request.
- ☒ Complete the Service Engineer Comments section if there are additional comments.
- ☒ Review this service, parts replaced, and test results obtained with the customer.
- ☒ If the instrument firmware was updated, record the details of the change in the Service Engineer's Comments box. Systems in a compliant environment may need additional documentation.
- ☒ **Complete the Signature Page with both Service Engineer and Customer signatures.**

Test Results

Instrument Performance Test Results Table

Note: These measurements do not form part of any specification and are for reference only.

	Pre PM Sensitivity Check		Post PM Sensitivity Check	
	Radial	Axial *	Radial	Axial*
Zn 213.857 nm SRBR	1783.8	2579.2	2230.8	3562.9
Mn 257.610 nm SRBR	9870.5	19614.9	10883.0	17884.4
Al 396.152 nm SBR	5.2	8.0	7.0	7.4
K 766.491 nm SBR	3.0	14.6	3.4	8.7

* Axial result is not applicable for G8016AA, G8012AA Radial View instruments.

Instrument Test Results Table

Note: The Instrument Test results are for systems using ICP Expert version 7.3 and above only.

Instrument Test	Result
Subsystem Communications Test	Pass
Air Flow	Pass
Water Flow	Pass
Gas Flows	Pass
RF Generator	Pass
Camera Test	Pass
Optics Test	Pass
Nebulizer test	Pass

ICP-OES Status Results Table

Note: These measurements do not form part of any specification and are for reference only.

Measurement	Standby Mode		Plasma On	
Mains Voltage	231.288	VAC	226.380	VAC
Mains Current	0.083	A	0.100	A
Instrument Temperature	21.1	°C	21.6	°C
RF Air Flow (sensor speed)	44.0	Hz	41.1	Hz
Plasma Exhaust Temperature	No measurement		47.1	°C
Water Flow Oscillator	No measurement		1.14	L/min
Water Flow Detector	0.92	L/min	0.90	L/min
Water Inlet Temperature	19.5	°C	18.5	°C
Polychromator Temperature	35.4	°C	36.4	°C
CCD Temperature	-40.1	°C	-39.9	°C
Thermal Stabilizer	35.0	°C	35.0	°C
Argon Supply Pressure	634.43	kPa	583.77	kPa
Purge Gas Supply Pressure*1	631.28	kPa	604.02	kPa
Option Gas Supply Pressure*1	-	kPa	-	kPa
Nebulizer Flow	No measurement		0.70	L/min
Nebulizer Back Pressure	No measurement		273.61	kPa
Plasma Gas Flow	No measurement		11.96	L/min
Auxiliary Gas Flow	No measurement		1.00	L/min
RF Power	No measurement		1200.4	W
RF Supply Current	No measurement		8.228	A
RF Supply Voltage	No measurement		194.529	V

*1 If option installed

Consumed PM Parts

Part Description	Part Number	Product or Model# where used	Quantity consumed
Axial Pre-Optic Window	G8010-68014	G8010A, G8011A, G8014A/G8015A	1
Radial Pre-Optic Window	G8010-68015	All	1
Agilent Cool Clear Coolant Fluid	5799-0037	Agilent Water Recirculator	1
Purge Gas Filter	G8010-60136	All	1
Air inlet filter	G8000-68002	All	1
High Capacity Air Filter	G8010-60189	Optional	1
Rotor seal for 6-7 port valve for AVS6/7	G8494-60002	G8494A/G8495	1
Rotor seal for 4 port valve for AVS4	G8493-60002	G8493A	1
Rinse solution to rinse station 2.5mm id x 1m	G8410-80123	SPS 4	1
Barb connector 2.5mm-1.5mm ID	G8410-80124	SPS 4	1
PVC waste tubing, 8mm od x 5mm id, 2m	G8410-80122	SPS 4	1
Additional Parts may be required from engineer's stock:			
X axis drive belt	5410047500	SPS 3	1
Z axis drive belt	5410047400	SPS 3	1
Peristaltic pump tubing, PVC SolvaFlex, 3 bridged,	3710049000	SPS 4	1

Consumed Parts Reference
(Purchased by customer, not included as part of PM)

☒ Section Not Applicable.

Part Description	Part Number	Product or Model# where used	Quantity consumed
------------------	-------------	------------------------------	-------------------

Signature Page

Service Engineer Comments (optional)

If there are any specific points you wish to note as part of performing the installation or other items of interest for the customer, please write in this box.

Service Verification

Service Request Number:

6008 006610

Service Engineer Name:

Worawit T.

Service Engineer Signature:

Worawit T.

Total number of pages in this document:

14

Date Service Completed:

3 Nov 2025

Customer Name:

Aphorn O.


Customer Signature:

Aphorn O.


CERTIFICATE OF CALIBRATION

Certificate No. : SP25-021

Page 1 of 5

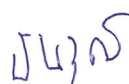Customer : United Analyst and Engineering Consultant Co.,Ltd. (Head Office)**Address :** 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260**Location of calibration :** Instrument room (207)**Equipment :** UV-Vis Spectrophotometer**Manufacturer :** HITACHI**Model :** U-2900**Serial No. :** 21E22-009**ID No. :** UAE.WAT.051/2564**Received Date :** 26 May 2025**Calibration Date :** 26 May 2025**Issue Date :** 29 May 2025**Condition Instrument :** Good**Calibrated by :**
(Mr.Tanawut Rittidach)

Technical Manager

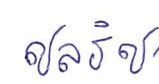Approved by :
(Ms.Chonticha Sangngern)

Quality Manager

The calibration result is applied only to the above calibrated item and was found accurate as shown on date and place of calibration only.

The measurement capability of the laboratory and its traceability to recognized national standards and to the unit of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the DQE Services Co., Ltd.

REPORT OF CALIBRATION

Certificate No. : SP25-021

Page 2 of 5

Environment Condition : Ambient Temperature 25 ± 5 °CRelative humidity 55 ± 20 %RH**Calibration method :** In-house method CP-01 Based on ASTM E275-08**Certified Reference Materials :**

Material	Serial No.	Certificate No.	Due date
Absorbance Standard set	25760	115663	25 October 2025
Absorbance Standard set	25757	115638	25 October 2025
Wavelength Standard set	25806	115657	25 October 2025
Wavelength Standard set	25758	115665	25 October 2025

Traceability : This certification is traceable to the International System of Unit maintained at National -

Institute of Standards and Technology (NIST) through Starna Scientific Limited

Spectral Band Width of UUC : 1.5 nm.**Scan Speed of UUC :** 200 nm/min**Scan Interval of UUC :** 0.1 nm.**Resolution of UUC :** Photometric 0.001 Abs.

Wavelength 0.1 nm.

REPORT OF CALIBRATION

Certificate No. : SP25-021

Page 3 of 5

Calibration Results : Without adjustment

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor <i>k</i>
420	0.0000	0.000	0.0000	0.0028	2.00
	0.5780	0.574	0.0040	0.0031	2.00
	1.0484	1.043	0.0054	0.0029	2.00
	2.1876	2.185	0.0026	0.0075	2.00
440	0.0000	0.000	0.0000	0.0028	2.00
	0.5595	0.558	0.0015	0.0034	2.00
	1.0239	1.022	0.0019	0.0035	2.00
	2.1230	2.121	0.0020	0.0079	2.00
465	0.0000	0.000	0.0000	0.0028	2.00
	0.5230	0.519	0.0040	0.0030	2.00
	0.9633	0.960	0.0033	0.0029	2.00
	1.9753	1.973	0.0023	0.0070	2.00
546.1	0.0000	0.000	0.0000	0.0028	2.00
	0.5181	0.516	0.0021	0.0031	2.00
	1.0002	0.998	0.0022	0.0033	2.00
	1.9973	1.995	0.0023	0.0084	2.00
590	0.0000	0.000	0.0000	0.0028	2.00
	0.5517	0.550	0.0017	0.0030	2.00
	1.0803	1.080	0.0003	0.0030	2.00
	2.0373	2.036	0.0013	0.0079	2.00
635	0.0000	0.000	0.0000	0.0028	2.00
	0.5591	0.558	0.0011	0.0031	2.00
	1.0518	1.051	0.0008	0.0030	2.00
	1.9274	1.926	0.0014	0.0079	2.00

เอกสารไม่ควบคุม

REPORT OF CALIBRATION

Certificate No. : SP25-021

Page 4 of 5

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor <i>k</i>
235	0.0000	0.000	0.0000	0.0050	2.00
	0.7469	0.748	-0.0011	0.0057	2.00
257	0.0000	0.000	0.0000	0.0050	2.00
	0.8674	0.866	0.0014	0.0059	2.00
313	0.0000	0.000	0.0000	0.0050	2.00
	0.2919	0.291	0.0009	0.0051	2.00
350	0.0000	0.000	0.0000	0.0050	2.00
	0.6430	0.640	0.0030	0.0055	2.00

REPORT OF CALIBRATION

Certificate No. : SP25-021

Page 5 of 5

Wavelength Accuracy :

CRMs Values (nm.)	UUC Reading (nm.)	Correction (nm.)	Uncertainty (nm.)	Coverage factor <i>k</i>
241.72	241.0	0.72	0.18	2.00
279.45	278.8	0.65	0.18	2.00
287.81	287.2	0.61	0.18	2.00
334.06	333.6	0.46	0.18	2.00
360.93	360.6	0.33	0.18	2.00
418.59	418.2	0.39	0.18	2.00
445.94	445.6	0.34	0.18	2.00
460.02	459.8	0.22	0.18	2.00
536.59	536.6	-0.01	0.18	2.00
637.98	638.0	-0.02	0.18	2.00
431.38	431.2	0.18	0.18	2.00
472.50	472.4	0.10	0.18	2.00
513.47	513.4	0.07	0.18	2.00
528.88	528.8	0.08	0.18	2.00
573.17	573.2	-0.03	0.18	2.00
585.35	585.0	0.35	0.20	2.00
684.40	684.4	0.00	0.18	2.00
740.72	741.0	-0.28	0.20	2.00
748.55	748.8	-0.25	0.18	2.00
807.03	807.2	-0.17	0.18	2.00
879.28	879.6	-0.32	0.18	2.00

Remark : - UUC = Unit Under Calibration

- N/A = Not Available

- The result expanded uncertainty of measurement U is stated as the standard uncertainty of measurement multiplied by the coverage factor k ,

which for a normal distribution corresponds to a coverage probability of approximately 95%

- End of Certificate -

เอกสารไม่ควบคุม



Certificate of Calibration

Cert.No.: 26CH144

Page.: 1 of 3

Equipment : Conductivity Meter
Manufacturer : Horiba
Model : LAQUA-EC210
Serial No. : HC9L0013
ID No. : UAE.EFM.011/2563(EFM.SCT.05/63)
Condition As-Received: Used Item
Received Date : 02 February 2026
Calibration Date : 03 February 2026
Reference : 2602-0038WSC-3
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road, Bangchak,
Phrakhanong, Bangkok 10260
Ambient Temperature : $(25 \pm 2.5) ^\circ\text{C}$
Relative Humidity : $(50 \pm 15) \%$
Calibration Procedure: In -house method :
- CP-CH6 by direct measurement
with certified reference material (CRM)
- CP-CH8 by comparison with temperature standard
Calibrated by : Walalak Sirithean
Approved by : _____
Approved Signatory
() Chakrit Waewwanjua
() Ponpan Paipim
(✓) Saithip Meangmai
Issue Date : 04 February 2026

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Cert.No.: 26CH144

Page.: 2 of 3

Condition of this result of calibration

1. Reference Standard Instrument :-

<u>Instrument</u>	<u>Serial No.</u>	<u>ID No.</u>	<u>Certificate No.</u>	<u>Due date</u>
1) Thermometer	9549224	130RC003	25I440	16 Apr 2026
2) Ref. Std.Thermometer	4982054	110RC044	25I708	03 July 2026

- This measurement result is traceable to SI through Technology Promotion Association (Thailand - Japan)

2. Certified Reference Materials :-

- Conductivity calibration solution, CPA chem Ltd., The measurement results are traceable to SI through CPA chem Ltd., ANSI-ASQ National Accreditation Board, Accredited No. AR-1835

<u>Conductivity Solution</u>	<u>Manufacturer</u>	<u>Lot No.</u>	<u>Exp. date</u>
1413.5 μ S/cm	CPA Chem	1150497	30 Oct 2026
12.881 mS/cm	CPA Chem	1081099	11 Mar 2026

- Control Conductivity calibration solution temperature by Water bath (25 ± 0.1) °C

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration results

Function : Conductivity Measurement

(*) After Adjustment at 1413.5 μ S/cm

Conductivity Electrode Serial No.: 9B9F0286

Standard Conductivity Solution	Before Adjustment UUC* Reading	After Adjustment UUC* Reading	Uncertainty of Measurement ($\pm$)	Coverage factor <i>k</i>
1413.5 μ S/cm	1431 μ S/cm	1413 μ S/cm	9.2 μ S/cm	2.00
12.881 mS/cm	12.92 mS/cm	12.76 mS/cm	0.086 mS/cm	2.00

Remark : - UUC* = Unit Under Calibration



Cert.No.: 26CH144

Page.: 3 of 3

Calibration Results

Function : Temperature Measurement

(*) Without adjustment

This equipment was connected with Temperature Probe;

- Model : 9383
- Serial No. : 9B9F0286

Dimension of probe;

- Length : 110 mm
- Diameter : 16 mm
- Immersion Depth : 100 mm

Calibration Point (°C)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty of Measurement (± °C)	Coverage factor <i>k</i>
15.0	15.002	15.0	-0.002	0.13	2.00
30.0	30.001	30.0	-0.001	0.13	2.00
45.0	45.001	45.0	-0.001	0.13	2.00

Remark : - UUC* = Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3 : EQUIPMENT CALIBRATION AND TESTING SERVICES

534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250

TEL. 0-2717-3000 FAX. 0-2719-9484

Certificate of Testing

Cert.No.: 25TW213

Page.: 1 of 2

Equipment :	DO Meter
Manufacturer :	YSI
Model :	4010-2W
Serial No. :	20260326
ID No. :	UAE.WAO.060/2563
Received Date :	17 October 2025
Test Date :	22 October 2025
Reference :	2510-0505DSC-1
Submitted by :	United Analyst and Engineering Consultant Co.,Ltd. 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260
Laboratory Condition :	Temperature (25 ± 5) °C Humidity (50 ± 20) %
Test Procedure :	In - house method : CP-CH9 by Comparison Technique with Azide Modification Method
Tested by :	Walalak Sirithean 
Approved by :	Approved Signatory
() Chakrit Waewwanjua () Ponpan Paipim (✓) Saithip Meangmai	
Issue Date :	23 October 2025

เอกสารไม่ควบคุม



Cert.No.: 25TW213

Page.: 2 of 2

Condition of this result of calibration

1. Reference Standard Instruments :

This measurement result is traceable to the International System of Unit through the reference standards laboratory of Industrial Calibration Center, Technology Promotion Association (Thailand-Japan).

<u>Instruments</u>	<u>ID No.</u>	<u>Certificate No.</u>	<u>Due Date</u>
1. Burette	130BU10	25CG1126	18 Mar 2027
2. Balance	110RC001	25MM316	02 July 2026

2. Standard Material :-

<u>Material</u>	<u>Manufacturer</u>	<u>Lot.No.</u>	<u>Assay</u>
Sodium Thiosulfate 5-Hydrate AR	KEMAUS	2203162447	99.6%

Result : Dissolved Oxygen Meter Adjustment With Air 100 %

Dissolved Oxygen Probe No.: 22M102385

Titration Method (Azide Modification Method) (mg/L)	DO Meter Reading (mg/L)	Standard Deviation (mg/L)
8.22	8.31	0.0055

This report was certified only for the instrument we tested. It is allowable to use for study
Intend to use for advertising and referral purpose is prohibited. This report may not be reproduced
other in full, without written approval of the laboratory

-o0o-

เอกสารไม่ควบคุม

Calibration Certificate

Certificate No.: 2601956-001-01

Page : 1 of 2

Equipment: pH Meter
Manufacturer: METTLER TOLEDO
Model: SevenEasy pH
Serial No.: 1231155210
ID No.: UAE.WAT.010/2553
Order No.: 2601956
Operation No.: 2601956-001
Date of Receipt: 20 February 2026
Date of Calibration: 23 February 2026

Customer name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.

Address: 3 Soi Udomsuk 41, Sukhumvit Road,
Bangchack, Prakhnong, Bangkok 10260

Metrological traceability:

1. Calibration Method W-CC-002 : In house method based on direct measurement by using standard voltage calibrator and certified reference material (CRM)
2. Reference Standards / Certified Reference Material

Instruments	Serial / ID No.	Manufacturer	Certificate No.	Due Date
2.1 DC Voltage Calibrator	2709007	Fluke	25E1753	22 May 2026
2.2 Digital Thermometer	2709007	Fluke	2600328-002-01	25 October 2026
2.3 Thermo-Hygro Meter	NFI.BTH 014/23	testo	2503232-014-01	10 June 2026
Certified Reference Material	Lot. No.	Manufacturer	Ref N	Expire Date
2.4 pH buffer 4.008 (Primary pH buffer Solution)	1016435	CPAchem	PH216.L5	25 July 2026
2.5 pH buffer 6.865 (Primary pH buffer Solution)	1016432	CPAchem	PH217.L5	25 July 2026
2.6 pH buffer 10.01 (Primary pH buffer Solution)	1088003	CPAchem	PH220.L5	18 March 2026
2.7 pH buffer 7.00 (Standard pH buffer Solution)	C03282	HACH LANGE GmbH	S11M004	11 June 2027

3. This certification is traceable to The International System of Unit (SI Unit)
 - 3.1 Instruments No.2.1 through NSC-ONSC Laboratory Accreditation of Calibration No.0008 (TPA)
 - 3.2 Instruments No.2.2 and 2.3 through NSC-ONSC Laboratory Accreditation of Calibration No.0061 (NFI)
 - 3.3 Certified Reference Material No.2.4 to 2.6 traceable to CPAchem Ltd is accredited to ISO 17034 (ANSI accredited No. AR-1835)
 - 3.4 Certified Reference Material No.2.7 traceable to Hach Lange GmbH is accredited to DIN EN ISO 17034 (DAkkS accredited No.D-RM-15184-01-00)
4. This certificate was certified only for the instrument we calibrated.
5. This result of calibration was found accurate as shown on date and place of calibration only.

Calibrated by Mr.Manas Somsak
Specialist

Approved by



(Mr.Pheraphat Tuanjit)

Manager, Division of Calibration Laboratory

Responsible for the Technical Management Team

Date of Issue: 2 March 2026

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 02 Date: 01-10-68



Calibration Certificate

Certificate No.: 2601956-001-01

Page : 2 of 2

Calibration Condition:

Place of Calibration : Chemical Calibration Laboratory, National Food Institute
Environment Condition: Ambient Temperature (24.0 ± 0.2) °C Relative Humidity: (49.8 ± 2) %
Calibration Item Condition: Used Item Type: Bench top
Range: 0 - 14 pH Resolution: 0.01 pH, 1 mV

Calibration results:

1. Calibration of pH Meter (Manual Temperature Compensation at 25 °C)

Nominal pH	DC Voltage Standard (mV)	Average Indicator Reading		Uncertainty ($\pm$ mV)	Coverage Factor (k)
		mV	pH		
0	414.119	414	-0.01	0.58	2.00
2	295.813	296	1.99	0.58	2.00
4	177.463	178	4.00	0.58	2.00
6	59.159	59	6.00	0.58	2.00
7	0.000	0	7.00	0.58	2.00
8	-59.159	-59	8.00	0.58	2.00
10	-177.462	-177	10.00	0.58	2.00
12	-295.813	-296	12.00	0.58	2.00
14	-414.119	-414	14.00	0.58	2.00

2. Calibration of pH Meter with Electrode (Manual Temperature Compensation at 25 °C)

Equipment: pH Electrode Manufacturer: METTLER TOLEDO Model: InLab Solids
Serial No.: C507554932 ID.No. N/A Type: Combined Electrode

Performance of Electrode system (Three-Point Calibration at pH 4, 7 and 9)

Certified Value @25 °C (pH)	Average Indicator Reading		Relative Slope (%)	Uncertainty ($\pm$ pH)	Coverage Factor (k)
	pH	mV			
4.008	4.01	184	-	0.0071	2.00
7.000	7.00	8	99.5	0.0086	2.00
10.010	10.01	-166	97.7	0.0085	2.00
6.876	6.88	15	-	0.0071	2.00

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k= 2, providing a level of confidence of approximately 95 %.

----- End of Report -----

F-CS-009 Revision: 02 Date: 01-10-68



Calibration Certificate

Certificate No.: 2601956-002-01

Page : 1 of 2

Equipment: Digital Thermometer with RTD (pH Meter)
Manufacturer: METTLER TOLEDO **Customer Name:** UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Model: SevenEasy pH **Address:** 3 Soi Udomsuk 41, Sukhumvit Road,
Serial No.: 1231155210 Bangchack, Prakanong, Bangkok 10260
ID No.: UAE.WAT.010/2553
Order No.: 2601956
Operation No.: 2601956-002
Date of Receipt: 20 February 2026
Date of Calibration: 23 February 2026

Metrological traceability:

1. Calibration Method:

- NFI Method W-TE-025 : based on comparison with a standard thermometer.
- TM-CG01-TT01-V04: CALIBRATION OF DIGITAL THERMOMETER WITH TEMPERATURE SENSORS.
- The temperature scale used in this laboratory is the International Temperature Scale of 1990 (ITS-90).

2. Reference Standard Instrument:

Instrument	Model / Type	Serial No. / ID No.	Certificate No.	Due Date	Through
HANDHELD THERMOMETER	1523	2118154	PSL-T 0926/68	26-Jun-26	TISTR
Platinum Resistance Thermometer (PRT)	5627A	877332			

Support Equipment : - Low Temperature Bath (AMETEK RTC-187), Model: RTC-187C , S/N: 670930-00018

- This certificate is traceable to the International System of Units (SI).
- This certificate is valid only for the instrument that was calibrated.
- The result of this calibration was found to be accurate only on the date and at the place of calibration shown.

Calibrated by Mr.Manas Somsak
Specialist

Approved by 
(Mr.Pheraphat Tuanjit)

Date of Issue: 2 March 2026

Manager, Division of Calibration Laboratory
Responsible for the Technical Management Team

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 02 Date: 01-10-68



Calibration Certificate

Certificate No.: 2601956-002-01

Page : 2 of 2

Calibration Condition :

Equipment: Digital Thermometer with RTD
UUC* Resolution: 0.1 °C
Condition of the item: Normal & Good
Without adjustment / AS FOUND

Place of Calibration: Chemical Calibration Laboratory

National Food Institute

Environment Conditions: Ambient Temperature (23.1 ± 0.3) °C
Relative Humidity (49.8 ± 2) %
Line Voltage (220 ± 10) V.

Calibration results:

- Temperature Measurement Results

Calibration Point (°C)	Standard Thermometer Value (°C)	UUC* Reading Value (°C)	UUC* Correction Value (°C)	Measurement Uncertainty (± °C)
20	20.002	20.1	-0.1	0.098
25	25.001	25.1	-0.1	0.098
30	30.003	30.1	-0.1	0.098

Note

- UUC* : Unit Under Calibration
- The probe was immersed in a liquid bath or a dry bath to a minimum depth of 120 mm.
- Description of probe, model : N/A S/N : N/A
- Dimension of probe : Diameter 4 mm. Length 120 mm. Sheath material: Plastic

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor $k = 2$, providing a level of confidence of approximately 95 %

----- End of Report -----

F-CS-009 Revision: 02 Date: 01-10-68

Certificate of Calibration

Certificate No.: 260311-2-BL001-26

Code No.: BL001-26

Page: 1 of 4

Customer Name: United Analyst and Engineering Consultant Co., Ltd.
Address: 3 Soi Udom suk 41, Sukhumvit Rd., Bang Chak, Phar Khanong, Bangkok 10260

Equipment: Electronic Balance
Manufacturer: Mettler Toledo
Model: XSR205DU
Serial No.: C210685394
Asset No. : UAE.WAO.010/2565

Building : UAE number 3 Floor : 2 Room : 208

Received Date: March 11, 2026

Date of Calibration : March 12, 2026

Calibration Conditions: Temperature 23.0 °C to 23.6 °C
Humidity 51.5 % to 54.4 %
Pressure 755.7 mmHg to 757.9 mmHg

Calibrated by: Thanachat Chatpahol

Approved by: Suwit Chotnok

Signature: 

Issued Date: March 16, 2026

Note : 1) The Uncertainties are for a confidence probability of approximately 95%

2) This Certificate is valid only to the item calibrated on date and place of calibration.

3) This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the United Analyst and Engineering Consultant Co.,Ltd. (UAE)

เอกสารไม่ควบคุม

Certificate No.: 260311-2-BL001-26

Code No.: BL001-26

Page: 2 of 4

Equipment: Electronic Balance

Manufacturer: Mettler Toledo

Model: XSR205DU

Readability: 0 to 81 (g): 0.01 mg , > 80 to 200 (g) : 0.1 mg

Serial No.: C210685394

ID No.: UAE.WAO.010/2565

Max. Capacity: 81 g / 220 g

Calibration Date: March 12, 2026

Condition As-Received: In Condition

Condition of Equipment:

Condition of This Result of Calibration:

1. Calibration Method: This instrument was calibrated by method UAE.CP.CAL.006 In-House Method based on UKAS Lab 14 : 2022

2. Reference Standards:

Reference Standard:	Model	Serial No.	Calibrated By	Certificate No.	Traceability	Due Date
Standard Weight Class E2 (OIML)	1 mg to 1 kg	C510641605	METTLER TOLEDO GmbH	C510641605	METTLER TOLEDO GmbH	3-Mar-27
Standard Weight Class E2 (OIML)	1 mg to 1 kg	B749109122	AMARC	25-009359	AMARC	21-Jan-27

Instrument	Model	Serial No.	Calibrated By	Certificate No.	Traceability	Due Date
Thermo-Hygro-Baro Meter	MHB-382SD	AK.46457	SUCCESS	SG-H-00992/68	Success Gateway	26-Nov-26
Thermo-Hygro-Baro Meter	MHB-382SD	AK.46457	TPA	25P4464	TPA	18-Dec-26

3. This certification is traceable to SI Unit

4. This certification was certified only for the instrument we calibrated

5. This result of calibration was found accurate as shown on date and place of calibration only.

6. Through the reference standard laboratory of Mettler-Toledo GmbH C510641605, SAS accreditation number SCS 0032

7. Through the reference standard laboratory of AMARC 25-009359, TISI accreditation number 0152

Calibration Result:




1. Repeatability of Reading:

Nominal Value (g)	Standard Deviation of Reading (g)
80	0.000019
200*	0.000055

2. Eccentric or off-center loading

A mass of 100 g was placed and moved to various positions on pan

The Balance reading obtained is given in the table.

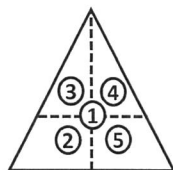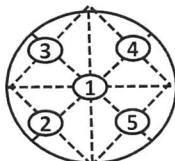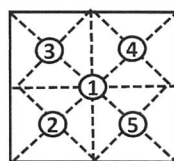

1 (g)	2 (g)	3 (g)	4 (g)	5 (g)	Maximum Difference (g)
100.0002	100.0002	100.0003	100.0004	100.0003	0.0001

e.aj

เอกสารไม่ควบคุม

Certificate No.: 260311-2-BL001-26

Code No.: BL001-26

Page: 3 of 4

Equipment:	Electronic Balance	Manufacturer:	Mettler Toledo
Model:	XSR205DU	Readability:	0 to 81 (g): 0.01 mg , > 80 to 200 (g) : 0.1 mg
Serial No.:	C210685394	ID No.:	UAE.WAO.010/2565
Max. Capacity:	81 g / 220 g		
Calibration Date:	March 12, 2026		

Calibraton Result: (Continued)

Calibraton Range: 0 to 80 (g)

Calibraton Adjustment: Internal Calibration

3. Error of indication from nominal or conventional mass value: (Range: 0 to 81 (g) ; Resolution: 0.00001 g)

Nominal Value (g)	Reference Value (g)	Indication (g)	Correction (g)	Uncertainty (± mg)	Coverage Factor k
Unload	0.0000000	0.00000	0.00000	0.032	2.65
0.001	0.0010020	0.00101	0.00000	0.032	2.65
0.005	0.0050025	0.00501	-0.00001	0.032	2.65
0.01	0.0100005	0.01001	-0.00001	0.032	2.65
0.05	0.0500003	0.05001	-0.00001	0.032	2.65
0.10	0.1000083	0.10000	0.00001	0.033	2.65
0.50	0.5000121	0.50001	0.00000	0.033	2.65
1	1.0000119	1.00003	-0.00002	0.032	2.52
2	2.0000125	2.00003	-0.00002	0.033	2.52
5	5.0000146	5.00003	-0.00002	0.033	2.32
10	10.0000330	10.00006	-0.00003	0.034	2.21
20	20.000028	20.00006	-0.00003	0.041	2.07
30	30.000061	30.00017	-0.00011	0.055	2.00
40	40.000054	40.00021	-0.00016	0.066	2.00
50	50.000037	50.00023	-0.00019	0.068	2.00
60	60.000070	60.00029	-0.00022	0.084	2.00
80	80.000098	80.00039	-0.00030	0.12	2.00

e. cut

เอกสารไม่ควบคุม

Certificate No.: 260311-2-BL001-26

Code No.: BL001-26

Page: 4 of 4

Equipment:	Electronic Balance	Manufacturer:	Mettler Toledo
Model:	XSR205DU	Readability:	0 to 81 (g): 0.01 mg , > 80 to 200 (g) : 0.1 mg
Serial No.:	C210685394	ID No.:	UAE.WAO.010/2565
Max. Capacity:	81 g / 220 g		
Calibration Date:	March 12, 2026		

Calibration Result: (Continued)

Calibration Range: > 81 to 200 (g)

Calibration Adjustment: Internal Calibration

3. Error of indication from nominal or conventional mass value: (Range: > 81 - 220 (g) ; Resolution: 0.0001 g)

Nominal Value (g)	Reference Value (g)	Indication (g)	Correction (g)	Uncertainty ($\pm$ mg)	Coverage Factor k
Unload	0.0000000	0.0000	0.0000	0.032	2.65
90	90.0000910	90.0004	-0.0003	0.16	2.00
100	100.0000420	100.0003	-0.0003	0.16	2.00
110	110.0000750	110.0004	-0.0003	0.18	2.00
120	120.0000680	120.0004	-0.0003	0.19	2.00
130	130.0001010	130.0005	-0.0004	0.20	2.00
140	140.0000960	140.0004	-0.0003	0.21	2.00
150	150.0000790	150.0005	-0.0004	0.22	2.00
160	160.0001120	160.0004	-0.0003	0.23	2.00
170	170.0001050	170.0003	-0.0002	0.24	2.00
200	200.0001460	200.0001	0.0000	0.26	2.00

4. Effect of Tare test:

Tare Load (g)	Test Load (g)	Indication (g)	Correction (g)
100	20.000028	20.0001	0.0000
	40.000054	40.0001	0.0000
	60.000070	60.0000	0.0001
	80.000098	80.0000	0.0001
	100.000109	100.0000	0.0001

Remark: -

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k, providing a level of confidence of approximately 95%

o---o-End-o---o

e.c

เอกสารไม่ควบคุม

Certificate of Calibration

Certificate No.: 260318-1-BL007-26

Code No.: BL007-26

Page: 1 of 3

Customer Name: United Analyst and Engineering Consultant Co., Ltd.
Address: 3 Soi Udom suk 41, Sukhumvit Rd., Bang Chak, Phar Khanong, Bangkok 10260
Equipment: Electronic Balance
Manufacturer: Mettler Toledo
Model: AB204-S/FACT
Serial No.: 1129361010
Asset No. : UAE.WAS.002/2552
Building : UAE No.3 Floor : 1 Room : 107

Received Date: April 2, 2026

Date of Calibration : April 3, 2026

Calibration Conditions:	Temperature	start	21.9 °C	end	22.1 °C
	Humidity	start	60.6 %	end	60.1 %
	Pressure	start	756.6 mmHg	end	755.1 mmHg

Calibrated by: Thanachat Chatpahol

Approved by: Suwit Chotnok

Signature: 

Issued Date: April 7, 2026

Note : 1) The Uncertainties are for a confidence probability of approximately 95%
2) This Certificate is valid only to the item calibrated on date and place of calibration.
3) This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the United Analyst and Engineering Consultant Co.,Ltd. (UAE)

เอกสารไม่ควบคุม

Certificate No.: 260318-1-BL007-26

Code No.: BL007-26

Page: 2 of 3

Equipment: Electronic Balance

Manufacturer: Mettler Toledo

Model: AB204-S/FACT

Readability: 0.0001 g

Serial No.: 1129361010

ID No.: UAE.WAS.002/2552

Max. Capacity: 220 g

Calibration Date: April 3, 2026

Condition As-Received: In Condition

Condition of Equipment:

Condition of This Result of Calibration:

1. Calibration Method: This instrument was calibrated by method UAE.CP.CAL.006

In-House Method based on UKAS Lab 14 : 2022

2. Reference Standards:

Reference Standard:	Model	Serial No.	Calibrated By	Certificate No.	Traceability	Due Date
Standard Weight Class E2 (OIML)	1 mg to 1 kg	C510641605	METTLER TOLEDO GmbH	C510641605	Mettler-Toledo	3-Mar-27

Instrument	Model	Serial No.	Calibrated By	Certificate No.	Traceability	Due Date
Thermo-Hygro-Baro Meter	MHB-382SD	AK.46457	SUCCESS	SG-H-00992/68	Success Gateway	26-Nov-26
Thermo-Hygro-Baro Meter	MHB-382SD	AK.46457	TPA	25P4464	TPA	18-Dec-26

3. This certification is traceable to SI Unit

4. This certification was certified only for the instrument we calibrated

5. This result of calibration was found accurate as shown on date and place of calibration only.

6. Through the reference standard laboratory of Mettler-Toledo GmbH C510641605, SAS accreditation number SCS 0032

Calibration Result:




1. Repeatability of Reading:

Nominal Value (g)	Standard Deviation of Reading (g)
200*	0.00010

2. Eccentric or off-center loading

A mass of 100 g was placed and moved to various positions on pan

The Balance reading obtained is given in the table.

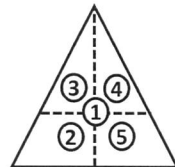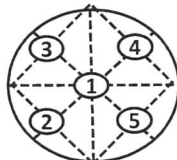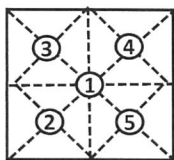

1 (g)	2 (g)	3 (g)	4 (g)	5 (g)	Maximum Difference (g)
100.0025	100.0022	100.0022	100.0029	100.0030	0.0005

ey

เอกสารไม่ควบคุม

Certificate No.: 260318-1-BL007-26

Code No.: BL007-26

Page: 3 of 3

Equipment:	Electronic Balance	Manufacturer:	Mettler Toledo
Model:	AB204-S/FACT	Readability:	0.0001 g
Serial No.:	1129361010	ID No.:	UAE.WAS.002/2552
Max. Capacity:	220 g		
Calibration Date:	April 3, 2026		

Calibraton Result: (Continued)

Calibraton Range: 0 - 200 g

Calibraton Adjustment: Internal Calibration

3. Error of indication from nominal or conventional mass value:

Nominal Value (g)	Reference Value (g)	Indication (g)	Correction (g)	Uncertainty ($\pm$ mg)	Coverage Factor k
Unload	0.0000000	0.0000	0.0000	0.17	2.37
0.01	0.0100005	0.0099	0.0001	0.17	2.32
0.05	0.0500003	0.0499	0.0001	0.17	2.32
0.1	0.1000083	0.0999	0.0001	0.17	2.32
0.5	0.5000121	0.5000	0.0000	0.17	2.32
1.0	1.0000119	1.0000	0.0000	0.17	2.32
10.0	10.000033	10.0003	-0.0003	0.17	2.32
40.0	40.000054	40.0010	-0.0010	0.18	2.23
50.0	50.000037	50.0013	-0.0012	0.18	2.23
80.0	80.000098	80.0020	-0.0019	0.19	2.11
100.0	100.000042	100.0025	-0.0024	0.20	2.10
120.0	120.000068	120.0028	-0.0027	0.22	2.06
150.0	150.000079	150.0035	-0.0035	0.24	2.04
160.0	160.000112	160.0036	-0.0035	0.25	2.03
200.0	200.000146	200.0049	-0.0047	0.28	2.00

4. Effect of Tare test:

Tare Load (g)	Test Load (g)	Indication (g)	Correction (g)
100	20.000028	20.0004	-0.0003
	40.000054	40.0009	-0.0008
	60.000070	60.0013	-0.0012
	80.000098	80.0017	-0.0016
	100.000042	100.0021	-0.0021

Remark: -

The report uncertainty of measurment was based on standard uncertainty multiplied by coverage factor k, providing a level of confidence of approximately 95%

o---o-End-o---o

เอกสารไม่ควบคุม

Certificate No. : C26030163

Page : 1 of 2

CERTIFICATE OF CALIBRATION

Equipment	: COD Test Tube Heater		
Meter Model	: HI839800-02	Serial No.	: 06480019101
Tube Heater Capacity	: 25 Vial Capacity	Resolution	: 0.1 °C
Temperature Range	: (-10 to 160)°C	Temperature of Reaction	: 150°C
Manufacturer	: Hanna Instruments	Made in	: Romania
Condition As-Received	: Used Product	Reference	: RE260497
Ambient Temperature	: (25 ± 2)°C	Relative Humidity	: (50 ± 15)% RH
Customer name	: United Analyst and Engineering Consultant Co., Ltd. 3 Soi Udomsuk 41, Sukhumvit Rd., Bangchak Phrakhanong, Bangkok 10260		
Received date	: 12 March 2026		
Calibrate date	: 17 March 2026		
Issue date	: 18 March 2026		
Calibrated Location	: Hanna Instruments (Thailand) Ltd.		
Calibration Procedure	: This calibrator was conducted by using in-house: calibration procedure CP-04 by using certified reference standard instruments.		

Calibrated by : Mr. Pichit Petthong**Approved by :**

(Mr. Anan Suwanchaisakul)

Authorized Signatory

Condition of this calibration result :

Reference standard instrument : This Certification is traceable to the SI Unit maintained through:

Instrument	Model	Serial No.	Certificate No.	Traceable
Digital Thermo-Hygrometer	HT-771SD	AI.07155	26H106	Technology Promotion Association (Thailand-Japan).
Data Acquisition Switch Unit	34970A	MY44065265	WK2507-171-1	WK Electric Co., Ltd.

Calibration Result :

Measurement Temperature Source Accuracy for COD Reactor.

Capacity	Nominal Value	Average Value	Error	Uncertainty ($\pm$)
25 Vials	150.0 °C	150.1 °C	0.1 °C	0.53 °C

Unit : °C

(1A)	(2A)	(3A)	(4A)	(5A)
149.7	149.7	150.1	150.4	150.4
(1B)	(2B)	(3B)	(4B)	(5B)
150.4	150.3	150.2	150.2	150.1
(1C)	(2C)	(3C)	(4C)	(5C)
149.9	150.2	150.0	150.2	149.9
(1D)	(2D)	(3D)	(4D)	(5D)
150.2	150.2	150.0	150.0	150.0
(1E)	(2E)	(3E)	(4E)	(5E)
150.1	149.8	150.2	150.2	149.8

Figure: Shows the location of the temperature source.

The report uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor $k=2$, providing a level of confidence of approximately 95%

** End of certificate **

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)

CORPORATE SERVICES 3 : EQUIPMENT CALIBRATION AND TESTING SERVICES

534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250

TEL. 0-2717-3000 FAX. 0-2719-9484

Certificate of Testing

Cert.No.: 26TW102

Page.: 1 of 2

Equipment :	DO Meter
Manufacturer :	Horiba
Model :	LAQUA-DO210
Serial No. :	HE9M0011
ID No. :	UAE.EFM.017/2563 (EFM.DO.06/63)
Received Date :	29 April 2026
Test Date :	30 April 2026
Reference :	2604-0723WSC-3
Submitted by :	United Analyst and Engineering Consultant Co.,Ltd. 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260
Laboratory Condition :	Temperature (25 ± 5) °C Humidity (50 ± 20) %
Test Procedure :	In - house method : CP-CH9 by Comparison Technique with Azide Modification Method
Tested by :	Walalak Sirithean 
Approved by :	Approved Signatory
() Chakrit Waewwanjua () Ponpan Paipim (✓) Saithip Meangmai	
Issue Date :	30 April 2026



Cert.No.: 26TW102

Page.: 2 of 2

Condition of this result of calibration

1. Reference Standard Instruments :

This measurement result is traceable to the International System of Unit through the reference standards laboratory of Industrial Calibration Center, Technology Promotion Association (Thailand-Japan).

<u>Instruments</u>	<u>ID No.</u>	<u>Certificate No.</u>	<u>Due Date</u>
1. Burette	130BU10	25CG1126	18 Mar 2027
2. Balance	110RC001	25MM316	02 July 2026

2. Standard Material :-

<u>Material</u>	<u>Manufacturer</u>	<u>Lot.No.</u>	<u>Assay</u>
Sodium Thiosulfate 5-Hydrate AR	KEMAUS	2203162447	99.6%

Result : **Dissolved Oxygen Meter Adjustment With Air 100 %**

Dissolved Oxygen Probe No.: 9K0H0035

Titration Method (Azide Modification Method) (mg/L)	DO Meter Reading (mg/L)	Standard Deviation (mg/L)
8.22	8.22	0.0089

This report was certified only for the instrument we tested. It is allowable to use for study
Intend to use for advertising and referral purpose is prohibited. This report may not be reproduced
other in full, without written approval of the laboratory

-o0o-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert. No.: 26LM81

Page.: 1 of 2

Equipment : DO Meter with Sensor

Manufacturer : Horiba

Model : LAQUA-DO210

Serial No. : HE9M0011

ID No. : UAE.EFM.017/2563 (EFM.DO.06/63)

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong, Bangkok 10260

Location : TPA On Site Calibration Laboratory

Received Order : 29 April 2026

Calibrated Date : 30 April 2026

Ambient Temperature : (26 ± 10) °C

Relative Humidity : (50 ± 30) %

AC Line Voltage : (220 ± 22) V

Calibrated by : Warakorn Lerngagtrakul

Approved by :

Kunchit

Approved Signatory

- () Chakrit Waewwanjua
() Suwit Imjai
(✓) Kunchit Promprat

Issue Date : 05 May 2026

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Equipment : DO Meter with Sensor
Condition As-Received : Used Item
Reference : 2604-0723WSC-4
Procedure Used :-

Cert. No.: 26LM81
Page.: 2 of 2

Calibration were conducted using in-house calibration procedure CP-OT01 according to comparison with Industrial Platinum Resistance Thermometer (IPRT) into Temperature Bath.

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

<u>Instrument</u>	<u>Serial No.</u>	<u>Cert. No.</u>	<u>Traceable</u>	<u>Due Date</u>
1) Digital Thermometer	2188080	25I953	TPA	05 Sep 2026

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This measurement result is traceable to the International System of Unit maintained through :

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function : Temperature measurement.

This instrument was connected with temperature sensor, S/N.: 9K0H0035

<u>Calibration Point</u> (°C)	<u>Immersion Depth</u> (mm)	<u>Standard Temperature</u> (°C)	<u>UUC* Reading</u> (°C)	<u>Error</u> (°C)	<u>Uncertainty</u> (± °C)	<u>Coverage Factor</u> <i>k</i>
15.0	80	15.004	15.0	-0.004	0.16	2.00
30.0	80	30.005	30.0	-0.005	0.16	2.00
45.0	80	45.004	44.9	-0.104	0.16	2.00

UUC* : Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-


CERTIFICATE OF CALIBRATION

Certificate No. : SP26-012

Page 1 of 5

Customer : United Analyst and Engineering Consultant Co.,Ltd. (Head Office)**Address :** 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260**Location of calibration :** Instrument room (207)**Equipment :** UV-Vis Spectrophotometer**Manufacturer :** HITACHI**Model :** U-2900**Serial No. :** 21E22-009**ID No. :** UAE.WAT.051/2564**Received Date :** 18 May 2026**Calibration Date :** 18 May 2026**Issue Date :** 21 May 2026**Condition Instrument :** Good**Calibrated by :**
(Mr.Pacit Mathavorn)

Calibration Engineer

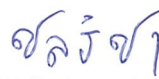Approved by :
(Ms.Chonticha Sangngern)

Quality Manager

The calibration result is applied only to the above calibrated item and was found accurate as shown on date and place of calibration only.

The measurement capability of the laboratory and its traceability to recognized national standards and to the unit of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the DQE Services Co., Ltd.

เอกสารไม่ควบคุม

FM-708-02 R01 1/11/2021

REPORT OF CALIBRATION

Certificate No. : SP26-012

Page 2 of 5

Environment Condition : Ambient Temperature 25 ± 5 °CRelative humidity 55 ± 20 %RH**Calibration method :** In-house method CP-01 Based on ASTM E275-08**Certified Reference Materials :**

Material	Serial No.	Certificate No.	Due date
Absorbance Standard set	25760	136858	5 November 2027
Absorbance Standard set	25757	136848	5 November 2027
Wavelength Standard set	25806	136824	5 November 2027
Wavelength Standard set	25758	136827	5 November 2027

Traceability : This certification is traceable to the International System of Unit maintained at National -

Institute of Standards and Technology (NIST) through Sarna Scientific Limited

Spectral Band Width of UUC : 1.5 nm.**Scan Speed of UUC :** 200 nm/min**Scan Interval of UUC :** 0.1 nm.**Resolution of UUC :** Photometric 0.001 Abs.

Wavelength 0.1 nm.

REPORT OF CALIBRATION

Certificate No. : SP26-012

Page 3 of 5

Calibration Results : Without adjustment

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor <i>k</i>
420	0.0000	0.000	0.0000	0.0028	2.00
	0.5772	0.575	0.0022	0.0031	2.00
	1.0478	1.043	0.0048	0.0029	2.00
	2.1879	2.197	-0.0091	0.0075	2.00
440	0.0000	0.000	0.0000	0.0028	2.00
	0.5581	0.558	0.0001	0.0034	2.00
	1.0224	1.022	0.0004	0.0035	2.00
	2.1227	2.131	-0.0083	0.0079	2.00
465	0.0000	0.000	0.0000	0.0028	2.00
	0.5225	0.518	0.0045	0.0029	2.00
	0.9628	0.961	0.0018	0.0028	2.00
	1.9744	1.981	-0.0066	0.0070	2.00
546.1	0.0000	0.000	0.0000	0.0028	2.00
	0.5174	0.515	0.0024	0.0031	2.00
	0.9990	0.997	0.0020	0.0033	2.00
	1.9953	1.995	0.0003	0.0084	2.00
590	0.0000	0.000	0.0000	0.0028	2.00
	0.5508	0.549	0.0018	0.0030	2.00
	1.0801	1.079	0.0011	0.0029	2.00
	2.0356	2.036	-0.0004	0.0079	2.00
635	0.0000	0.000	0.0000	0.0028	2.00
	0.5579	0.558	-0.0001	0.0031	2.00
	1.0512	1.050	0.0012	0.0029	2.00
	1.9257	1.928	-0.0023	0.0079	2.00

เอกสารไม่ควบคุม

REPORT OF CALIBRATION

Certificate No. : SP26-012

Page 4 of 5

Photometric Accuracy :

Wavelength (nm.)	CRMs Values (Abs)	UUC Reading (Abs)	Correction (Abs)	Uncertainty (Abs)	Coverage factor <i>k</i>
235	0.0000	0.000	0.0000	0.0050	2.00
	0.7455	0.742	0.0035	0.0056	2.00
257	0.0000	0.000	0.0000	0.0050	2.00
	0.8656	0.860	0.0056	0.0059	2.00
313	0.0000	0.000	0.0000	0.0050	2.00
	0.2910	0.288	0.0030	0.0051	2.00
350	0.0000	0.000	0.0000	0.0050	2.00
	0.6414	0.635	0.0064	0.0055	2.00

REPORT OF CALIBRATION

Certificate No. : SP26-012

Page 5 of 5

Wavelength Accuracy :

CRMs Values (nm.)	UUC Reading (nm.)	Correction (nm.)	Uncertainty (nm.)	Coverage factor <i>k</i>
241.72	241.3	0.42	0.18	2.00
279.45	279.0	0.45	0.18	2.00
287.81	287.4	0.41	0.18	2.00
334.06	333.8	0.26	0.18	2.00
360.93	360.5	0.43	0.18	2.00
418.59	418.0	0.59	0.18	2.00
445.94	445.5	0.44	0.18	2.00
453.66	453.5	0.16	0.18	2.00
460.02	459.5	0.52	0.18	2.00
536.59	536.5	0.09	0.18	2.00
637.98	638.0	-0.02	0.18	2.00
431.38	431.0	0.38	0.18	2.00
472.50	472.5	0.00	0.18	2.00
513.47	513.3	0.17	0.18	2.00
528.88	528.5	0.38	0.18	2.00
573.17	573.0	0.17	0.18	2.00
585.35	585.0	0.35	0.20	2.00
684.40	684.0	0.40	0.18	2.00
740.72	740.5	0.22	0.20	2.00
748.55	748.5	0.05	0.18	2.00
807.03	807.0	0.03	0.18	2.00
879.28	879.3	-0.02	0.18	2.00

Remark : - UUC = Unit Under Calibration

- N/A = Not Available

- The result expanded uncertainty of measurement U is stated as the standard uncertainty of measurement multiplied by the coverage factor k ,

which for a normal distribution corresponds to a coverage probability of approximately 95%

- * Indicates non TISI accredited

- End of Certificate -

เอกสารไม่ควบคุม

List of Instruments Certification for Water Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration	Remark
Water									
1	pH Meter	pH	Ecosence	pH100A 24H005156JEN	Technology Promotion Association (Thailand-Japan)	25CH1281	10 Nov 25	9 Nov 26	-
2	DO Meter	DO	Horiba	LAQUA-DO210 HE0H0003	Technology Promotion Association (Thailand-Japan)	25TW252	25 Dec 25	24 Dec 26	-



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert.No.: 25CH1281

Page.: 1 of 3

Equipment : pH Meter
Manufacturer : EcoSense
Model : pH100A
Serial No. : 24H005156JEN
ID No. : UAE.EFM.038/2567(EFM.pH.01/67)
Condition As-Received: Used Item
Received Date : 06 November 2025
Calibration Date : 10 November 2025
Reference : 2511-0154WSC-1
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road, Bangchak,
Phrakhanong, Bangkok 10260
Ambient Temperature : (25 ± 2.5) °C
Relative Humidity : (50 ± 15) %
Calibration Procedure : In - house method :
- CP-CH5 by direct measurement with DC voltage
standard and direct measurement with
certified reference material (CRM)
- CP-CH8 by comparison with temperature standard
Calibrated by : Walalak Sirithean

Approved by : 
Approved Signatory
() Chakrit Waewwanjua
() Ponpan Paipim
(✓) Saithip Meangmai
Issue Date : 11 November 2025

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Cert.No.: 25CH1281

Page.: 2 of 3

Condition of this calibration result

1. Reference Standard Instrument

<u>Instrument</u>	<u>Serial No.</u>	<u>ID No.</u>	<u>Cert. No.</u>	<u>Due Date</u>
1) Document Process Calibrator	54030049	130RC116	25E2743	25 Aug 2026
2) Ref. Standard Thermometer	4982054	110RC044	25I708	03 July 2026

- This measurement result is traceable to SI through Technology Promotion Association (Thailand - Japan)

2. Certified Reference Materials :The measurement results are traceable to SI through Hach Lenge GmbH Ltd.,
Deutsche Akkreditierungsstelle, Accredited No.D-RM-15184-01-00
:The measurement results are traceable to SI through CPA chem Ltd.,
ANSI-ASQ National Accreditation Board, Accredited No. AR-1835

<u>Buffer Solution</u>	<u>Manufacturer</u>	<u>Lot No.</u>	<u>Exp. date</u>
pH 4.007	CPA chem	1135353	16 Aug 2027
pH 6.999	Hach Lenge GmbH	C03273	13 May 2027
pH 10.010	CPA chem	1135355	16 Aug 2026

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration Results

Function : mV Measurement

Performing standard curve by Document Process Calibrator at pH (4,7)(7,10)

Unit Under Calibration	Nominal Value	Standard Voltage Input	Actual Reading		Uncertainty of Measurement (±mV)	Coverage factor <i>k</i>
	pH	mV	mV	pH		
pH Meter S/N.: 24H005156JEN	4.00	177.48	177	4.01	0.58	2.00
	7.00	0.00	0	7.00	0.58	2.00
	7.00	0.00	0	7.00	0.58	2.00
	10.00	-177.48	-177	10.01	0.58	2.00

เอกสารไม่ควบคุม



Cert.No.: 25CH1281

Page.: 3 of 3

Calibration Results

Function : pH Measurement

Performing three buffers standard curve by using buffer nominal pH (4,7)(7,10)

Unit Under Calibration	Standard pH Buffer Solution	Actual pH Reading	Actual mV Reading (mV)	Uncertainty of pH Measurement ($\pm$)	Coverage factor k
pH Electrode S/N.: 240904SIA605377	4.007	4.01	167	0.0079	2.00
	6.999	6.99	-5	0.0095	2.00
	6.999	7.00	-7	0.0092	2.00
	10.010	10.01	-180	0.0092	2.00

Function : Temperature Measurement

(*) Without adjustment

This equipment was connected with Temperature Probe;

- Model : -
- Serial No. : 240904SIA605377
- Dimension of probe
- Length : 110 mm.
- Diameter : 12 mm.
- Immersion Depth : 100 mm.

Calibration Point (°C)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty of measurement ($\pm$ °C)	Coverage factor k
15.0	15.002	15.0	-0.002	0.13	2.00
30.0	30.001	29.9	-0.101	0.13	2.00
45.0	45.002	44.7	-0.302	0.13	2.00

Remark - UUC* = Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k , providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3 : EQUIPMENT CALIBRATION AND TESTING SERVICES

534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250

TEL. 0-2717-3000 FAX. 0-2719-9484

Certificate of Testing

Cert.No.: 25TW252

Page.: 1 of 2

Equipment :	DO Meter
Manufacturer :	Horiba
Model :	LAQUA-DO210
Serial No. :	HE0H0003
ID No. :	UAE.EFM.083/2564(EFM.DO.02/64)
Received Date :	24 December 2025
Test Date :	25 December 2025
Reference :	2512-0641WSC-3
Submitted by :	United Analyst and Engineering Consultant Co.,Ltd. 3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260
Laboratory Condition :	Temperature (25 ± 5) °C Humidity (50 ± 20) %
Test Procedure :	In - house method : CP-CH9 by Comparison Technique with Azide Modification Method
Tested by :	Walalak Sirithean 
Approved by :	Approved Signatory
() Chakrit Waewwanjua	
() Ponpan Paipim	
(✓) Saithip Meangmai	
Issue Date :	29 December 2025

เอกสารไม่ควบคุม



Cert.No.: 25TW252

Page.: 2 of 2

Condition of this result of calibration

1. Reference Standard Instruments :

This measurement result is traceable to the International System of Unit through the reference standards laboratory of Industrial Calibration Center, Technology Promotion Association (Thailand-Japan).

<u>Instruments</u>	<u>ID No.</u>	<u>Certificate No.</u>	<u>Due Date</u>
1. Burette	130BU10	25CG1126	18 Mar 2027
2. Balance	110RC001	25MM316	02 July 2026

2. Standard Material :-

<u>Material</u>	<u>Manufacturer</u>	<u>Lot.No.</u>	<u>Assay</u>
Sodium Thiosulfate 5-Hydrate AR	KEMAUS	2203162447	99.6%

Result : Dissolved Oxygen Meter Adjustment With Air 100 %

Dissolved Oxygen Probe No.: 9K0E0260

Titration Method (Azide Modification Method) (mg/L)	DO Meter Reading (mg/L)	Standard Deviation (mg/L)
8.22	8.23	0.0045

This report was certified only for the instrument we tested. It is allowable to use for study
Intend to use for advertising and referral purpose is prohibited. This report may not be reproduced
other in full, without written approval of the laboratory

-o0o-

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL.0-2717-3000-29 FAX.0-2719-9484



Certificate of Calibration

Cert. No.: 26LM2

Page.: 1 of 2

Equipment : Do Meter with Sensor

Manufacturer : Horiba

Model : LAQUA-DO210

Serial No. : HE0H0003

ID No. : UAE.EFM.083/2564(EFM.DO.02/64)

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong, Bangkok 10260

Location : TPA On Site Calibration Laboratory

Received Order : 24 December 2025

Calibrated Date : 05 January 2026

Ambient Temperature : (26 ± 10) °C

Relative Humidity : (50 ± 30) %

AC Line Voltage : (220 ± 22) V

Calibrated by : Warakorn Lerngagtrakul

Approved by :

Kunchit

Approved Signatory

- () Chakrit Waewwanjua
() Suwit Imjai
(✓) Kunchit Promprat

Issue Date : 05 January 2026

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Equipment : Do Meter with Sensor
Condition As-Received : Used Item
Reference : 2512-0641WSC-4

Cert. No.: 26LM2
Page.: 2 of 2

Procedure Used :-

Calibration were conducted using in-house calibration procedure CP-OT01 according to comparison with Industrial Platinum Resistance Thermometer (IPRT) into Temperature Bath.

The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Digital Thermometer	2188080	25I953	TPA	05 Sep 2026

2. This certificate is valid only to the item calibrated on date and place of calibration.

3. This measurement result is traceable to the International System of Unit maintained through :

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment

Function : Temperature measurement.

This instrument was connected with temperature sensor, S/N.: 9K0E0260

Calibration Point (°C)	Immersion Depth (mm)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty (± °C)	Coverage Factor k
15.0	80	15.005	15.0	-0.005	0.16	2.00
30.0	80	30.006	30.0	-0.006	0.16	2.00
45.0	80	45.002	45.1	0.098	0.16	2.00

UUC* : Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor **k**, providing a level of confidence of approximately 95 %.

-o0o-