

ภาคผนวก ง
เอกสารสอบเทียบเครื่องมือ

List of Instrument Certificates for Environmental Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*
1	Analytical Balance	FAT OIL AND GREASE	Mettler Toledo	AB204-S/FACT / 1129361010	United Analyst and Engineering Consultant Co., Ltd.	260318-1-BL007-26	3/4/2026	2/4/2027
2	Analytical Balance	TOTAL DISSOLVED SOLIDS	Mettler Toledo	XSR205DU / C210685394	United Analyst and Engineering Consultant Co., Ltd (UAE)	260311-2-BL001-26	12/3/2026	11/3/2027
3	Analytical Balance	TOTAL SUSPENDED SOLIDS	Mettler Toledo	XSR205DU / C009071872	United Analyst and Engineering Consultant Co., Ltd. (UAE)	260311-2-BL002-26	12/3/2026	11/3/2027
4	DO Meter	BIOCHEMICAL OXYGEN DEMAND	YSI	5100 / 11B 101863	Technology Promotion Association (Thailand-Japan)	26TW50	2/3/2026	1/3/2027
5	Cooled Incubator	TOTAL COLIFORM BACTERIA	Binder	KB400 / WTB20200000015535	National Food Institute, Ministry of Industry, Thailand	2602303-006-01	11/3/2026	10/3/2027
6	pH Meter	pH	Horiba	LAQUA-PH210 / HA9M0047	technology promotion association (thailand-japan)	26CH231	26/2/2026	26/2/2027

Due Date of Calibration* : Based on the annual calibration plan. At least 1 time per year.



Certificate of Testing

Cert.No.: 26TW50
Page.: 1 of 2

Equipment : DO Meter
Manufacturer : YSI
Model : 5100
Serial No. : 11B 101863
ID No. : UAE.WAO.004/2554
Received Date : 27 February 2026
Test Date : 02 March 2026
Reference : 2602-0915DSC-1
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road, Bangchak,
Phrakhanong, Bangkok 10260
Laboratory Condition : Temperature (25 ± 5) °C
Humidity (50 ± 20) %
Test Procedure : In - house method : CP-CH9
by Comparison Technique with Azide Modification Method
Tested by : Walalak Sirithean
Approved by :
Approved Signatory
() Chakrit Waewwanjua
() Ponpan Paipim
(✓) Sathip Meangmai
Issue Date : 03 March 2026



Cert.No.: 26TW50
Page.: 2 of 2

Condition of this result of calibration

1. Reference Standard Instruments :
This measurement result is traceable to the International System of Unit through the reference standards laboratory of Industrial Calibration Center, Technology Promotion Association (Thailand-Japan).

Instruments	ID No.	Certificate No.	Due Date
1. Burette	130BU10	25CG1126	18 Mar 2027
2. Balance	110RC001	25MM316	02 July 2026

2. Standard Material :-

Material	Manufacturer	Lot.No.	Assay
Sodium Thiosulfate 5-Hydrate AR	KEMAUS	2203162447	99.6%

Result : Dissolved Oxygen Meter Adjustment With Air 100 %
Dissolved Oxygen Probe No.: 24F100202

Titration Method (Azide Modification Method) (mg/L)	DO Meter Reading (mg/L)	Standard Deviation (mg/L)
8.20	8.20	0.0084

This report was certified only for the instrument we tested. It is allowable to use for study
Intend to use for advertising and referral purpose is prohibited. This report may not be reproduced
other in full, without written approval of the laboratory

-oOo-

เอกสารไม่ควบคุม

เอกสารไม่ควบคุม



Cert.No.: 26CH231
Page.: 2 of 3

Certificate of Calibration

Cert.No.: 26CH231
Page.: 1 of 3

Equipment : pH Meter
Manufacturer : Horiba
Model : LAQUA-PH210
Serial No. : HA9M0047
ID No. : UAE.EFM.005/2563 (EFM.pH.05/63)
Condition As-Received: Used Item
Received Date : 24 February 2026
Calibration Date : 26 February 2026
Reference : 2602-0796WSC-2
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road, Bangchak,
Phrakhanong, Bangkok 10260
Ambient Temperature : (25 ± 2.5) °C
Relative Humidity : (50 ± 15) %
Calibration Procedure : In - house method :
- CP-CH5 by direct measurement with DC voltage
standard and direct measurement with
certified reference material (CRM)
- CP-CH8 by comparison with temperature standard
Calibrated by : Walalak Sirithean
Approved by :
Approved Signatory
() Chakrit Waewwanjua
() Ponpan Paipim
(✓) Sathip Meangmai
Issue Date : 27 February 2026

Condition of this calibration result

1. Reference Standard Instrument

Instrument	Serial No.	ID No.	Cert. No.	Due Date
1) Document Process Calibrator	54030049	130RC116	25E2743	25 Aug 2026
2) Ref. Standard Thermometer	4962054	110RC044	25I708	03 July 2026

- This measurement result is traceable to SI through Technology Promotion Association (Thailand - Japan)

2. Certified Reference Materials : The measurement results are traceable to SI through Hach Lenge GmbH Ltd.,
Deutsche Akkreditierungsstelle, Accredited No. D-RM-15184-01-00
The measurement results are traceable to SI through CPA chem Ltd.,
ANSI-ASQ National Accreditation Board, Accredited No. AR-1835

Buffer Solution	Manufacturer	Lot No.	Exp. date
pH 4.008	CPA chem	1187328	18 Jan 2028
pH 7.000	Hach Lenge GmbH	C03319	17 Oct 2027
pH 10.010	CPA chem	1187332	18 Jan 2027

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration Results

Function : mV Measurement

Performing standard curve by Document Process Calibrator at pH (4.7)(7.10)

Unit Under Calibration	Nominal Value	Standard Voltage Input	Actual Reading		Uncertainty of Measurement (±mV)	Coverage factor k
			mV	pH		
pH Meter S/N : HA9M0047	4.00	177.48	177.3	4.01	0.058	2.00
	7.00	0.00	0.0	7.02	0.058	2.00
	7.00	0.00	0.0	7.02	0.058	2.00
	10.00	-177.48	-177.2	10.01	0.058	2.00

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3 : Equipment Calibration and Testing Services.



Cert.No.: 26CH231
Page: 3 of 3

Calibration Results

Function : pH Measurement

Performing three buffers standard curve by using buffer nominal pH (4,7,10)

Unit Under Calibration	Standard pH Buffer Solution	Actual pH Reading	Actual mV Reading (mV)	Uncertainty of pH Measurement ($\pm$)	Coverage factor k
pH Electrode S/N: -	4.008	4.01	161.1	0.0079	2.00
	7.000	7.00	-10.4	0.0096	2.00
	7.000	7.00	-9.9	0.0095	2.00
	10.010	10.02	-177.9	0.0095	2.00

Function : Temperature Measurement

(*) Without adjustment

This equipment was connected with Temperature Probe;

- Model : -
- Serial No. : -
Dimension of probe
- Length : 103 mm.
- Diameter : 16 mm.
- Immersion Depth : 90 mm.

Calibration Point ($^{\circ}\text{C}$)	Standard Temperature ($^{\circ}\text{C}$)	UUC* Reading ($^{\circ}\text{C}$)	Error ($^{\circ}\text{C}$)	Uncertainty of measurement ($\pm$ $^{\circ}\text{C}$)	Coverage factor k
15.0	15.001	15.0	-0.001	0.13	2.00
30.0	30.001	30.0	-0.001	0.13	2.00
45.0	45.001	44.9	-0.101	0.13	2.00

Remark : UUC* = Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-o0o-



United Analyst and Engineering Consultant Co., Ltd.
3 Soi Udomsuk 41, Sukhumvit Road, Bangkok, Phra Khanong, Bangkok 10260
Tel 0 2763 2828 Fax 0 2763 2800 www.uaec consultant.com E-mail: uaec@uaec consultant.com



Certificate of Calibration

Certificate No.: 260311-2-BL002-26

Code No.: BL002-26

Page: 1 of 4

Customer Name: United Analyst and Engineering Consultant Co., Ltd.
Address: 3 Soi Udomsuk 41, Sukhumvit Rd., Bang Chak, Phra Khanong, Bangkok 10260.

Equipment: Electronic Balance
Manufacturer: Mettler Toledo
Model: XSR205DU
Serial No.: C009071872
Asset No.: UAE WAO 012/2563
Building : UAE number 3 Floor : 2 Room : 208

Received Date: March 11, 2026

Date of Calibration : March 12, 2026

Calibration Conditions: Temperature: 23.8 $^{\circ}\text{C}$ to 24.1 $^{\circ}\text{C}$
Humidity: 50.0 % to 59.0 %
Pressure: 755.5 mmHg to 755.6 mmHg

Calibrated by: Thanachat Chatpahol

Approved by: Suwit Chotnok

Signature:

Issued Date: March 16, 2026

Note : 1) The Uncertainties are for a confidence probability of approximately 95%

2) This Certificate is valid only to the item calibrated on date and place of calibration

3) This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the United Analyst and Engineering Consultant Co. Ltd. (UAE)

เอกสารไม่ควบคุม



United Analyst and Engineering Consultant Co., Ltd.
3 Soi Udomsuk 41, Sukhumvit Road, Bangkok, Phra Khanong, Bangkok 10260
Tel 0 2763 2828 Fax 0 2763 2800 www.uaec consultant.com E-mail: uaec@uaec consultant.com



Certificate No.: 260311-2-BL002-26

Code No.: BL002-26

Page: 2 of 4

Equipment: Electronic Balance
Model: XSR205DU
Serial No.: C009071872
Max. Capacity: 81 g / 220 g
Calibration Date: March 12, 2026
Condition As-Received: in Calibration

Manufacturer: Mettler Toledo
Readability: 0.1 to 81 g; 0.01 mg to 80 to 200 g; 0.01 mg
ID No.: UAE WAO 012/2563

Condition of Equipment:
Condition of This Result of Calibration:
1. Calibration Method: The instrument was calibrated by method UUC of CAL-006. In-house Method based on OIML R 110-2007

2. Reference Standards:
Reference Standard: Model Serial No. Calibrated By Certificate No. Traceability Due Date
Mettler Toledo Class 1 (FIM) 1 mg to 1 kg C510641605 METTLER TOLEDO GmbH 25-099359 Mettler Toledo 1-Mar-27
Mettler Toledo Class 1 (FIM) 1 mg to 1 kg B769109122 AMAR 25-099359 Mettler Toledo 27-Jan-27

Instrument: Model Serial No. Calibrated By Certificate No. Traceability Due Date
Thermo Hygro-Baro Meter M-19 38250 AR 46457 SUCC155 25-099359/68 Success Laboratory 26-Nov-26
Thermo Hygro-Baro Meter M-19 38250 AR 46457 T/A 25-099359 T/A 18-Dec-26

3. This certification is traceable to SI Unit.
4. This certification was certified only for the instrument we calibrated.
5. This result of calibration was found accurate as shown on date and place of calibration only.
6. Through the reference standard laboratory of Mettler Toledo GmbH C510641605. SAS accreditation number SC5 0032
7. Through the reference standard laboratory of AMAR 25-099359. TS accreditation number 0137

Calibration Results:


1. Repeatability of Reading:

Nominal Value (g)	Standard Deviation of Reading (g)
80	0.000008
200*	0.000045

2. Eccentric or off-center loading

A mass of 100 g was placed and moved to various position on pan.

The Balance reading obtained is given in the table:

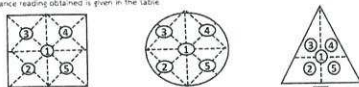

1	2	3	4	5	Maximum Difference (g)
(g)	(g)	(g)	(g)	(g)	
100.0000	100.0000	100.0000	100.0000	100.0000	0.0000

c. of

เอกสารไม่ควบคุม



United Analyst and Engineering Consultant Co., Ltd.
3 Soi Udomsuk 41, Sukhumvit Road, Bangkok, Phra Khanong, Bangkok 10260
Tel 0 2763 2828 Fax 0 2763 2800 www.uaec consultant.com E-mail: uaec@uaec consultant.com



Certificate No.: 260311-2-BL002-26

Code No.: BL002-26

Page: 3 of 4

Equipment: Electronic Balance
Model: XSR205DU
Serial No.: C009071872
Max. Capacity: 81 g / 220 g
Calibration Date: March 12, 2026

Manufacturer: Mettler Toledo
Readability: 0.1 to 81 g; 0.01 mg to 80 to 200 g; 0.01 mg
ID No.: UAE WAO 012/2563

Calibration Result: XSR205DU

Calibration Range: 0 to 801 g

Calibration Adjustment: Internal Calibration

3. Error of indication from nominal or conventional mass value: (Range: 0 to 81 g; Resolution: 0.00001 g)

Nominal Value (g)	Reference Value (g)	Indication (g)	Correction (g)	Uncertainty ($\pm$ mg)	Coverage Factor k
Unload	0.0000000	0.0000	0.0000	0.015	2.25
0.001	0.0010020	0.0010	0.0000	0.015	2.23
0.005	0.0050025	0.0050	-0.0002	0.015	2.23
0.01	0.0100005	0.0100	-0.0002	0.015	2.23
0.05	0.0500003	0.0500	-0.0001	0.015	2.21
0.10	0.1000008	0.1000	-0.0001	0.015	2.20
0.50	0.5000121	0.5000	-0.0002	0.016	2.15
1	1.0000119	1.0000	0.0000	0.016	2.12
2	2.0000125	2.0000	-0.0001	0.017	2.09
5	5.0000146	5.0000	-0.0002	0.020	2.05
10	10.0000330	10.0001	-0.0002	0.023	2.00
20	20.000028	20.0001	-0.0005	0.033	2.00
30	30.000061	30.0001	-0.0004	0.051	2.00
40	40.000054	40.0001	-0.0004	0.063	2.00
50	50.000037	50.0001	-0.0005	0.064	2.00
60	60.000070	60.0001	-0.0004	0.081	2.00
80	80.000098	80.0001	-0.0001	0.11	2.00

c. of

เอกสารไม่ควบคุม

Certificate No: 260311-2-BL001-26
Code No: BL001-26
Page: 4 of 4

Equipment: Electronic Balance
Model: XSR205DU
Serial No: C00901872
Max. Capacity: 81 g / 220 g
Calibration Date: March 12, 2026

Manufacturer: Mettler Toledo
Readability: 0.1 to 81 (g) 0.01 mg to 80 to 200 (g) 0.1 mg
ID No: UAE WAO 010/2563

Calibration Result: (Continued)
Calibration Range: > 80 to 200 (g)
Calibration Adjustment: Internal Calibration

3. Error of Indication from nominal or conventional mass value: (Range: > 80 to 200 (g) Resolution: 0.0001 g)

Nominal Value (g)	Reference Value (g)	Indication (g)	Correction (g)	Uncertainty ($\pm$ mg)	Coverage Factor k
Unload	0.0000000	0.0000	0.0000	0.02	2.25
90	90.0000910	90.0000	0.0000	0.16	2.00
100	100.0000420	100.0000	0.0000	0.16	2.00
110	110.0000750	110.0000	0.0000	0.17	2.00
120	120.0000680	120.0000	0.0001	0.18	2.00
130	130.0001010	130.0000	0.0000	0.20	2.00
140	140.0000960	140.0003	0.0001	0.21	2.00
150	150.0000790	150.0003	0.0001	0.21	2.00
160	160.0001120	160.0003	0.0001	0.23	2.00
170	170.0001050	170.0003	-0.0002	0.24	2.00
200	200.0001460	200.0004	-0.0002	0.26	2.00

4. Effect of Tare test

Tare Load (g)	Test Load (g)	Indication (g)	Correction (g)
100	20.000028	20.0000	0.0000
	40.000054	40.0002	-0.0001
	60.000070	60.0002	-0.0001
	80.000098	80.0002	-0.0001
	100.000042	100.0003	-0.0002

Remark:

The quoted uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k=2, providing a level of confidence of approximately 95%.

o--o-End-o--o

CU

เอกสารไม่ควบคุม

Calibration Certificate

Certificate No: 260333-006-01
Page: 2 of 2

Calibration Condition:

Equipment: CHAMBER (Incubator)
UUC* Resolution: 0.1 °C
Condition of the item: Normal & Good
Without adjustment: As Found
Time of record: 1 hour, 0 minute, 0 point

Place of Calibration: Microbiology Laboratory
UNITED ANALYST AND ENGINEERING CONSULTANT CO., LTD.
Environment Conditions: Ambient Temperature: $(17.2 \pm 1) ^\circ\text{C}$
Relative Humidity: $(51.5 \pm 1) \%$
Line Voltage: $(229 \pm 2) \text{ V}$

Calibration Results

Table 1: Reporting of Temperature

Calibration Point (°C)	Std.#1	Measured Temperature (°C) @ Sensor No. (Sensor No. 9 is REF)	Std.#2	Std.#3	Std.#4	Std.#5	Std.#6	Std.#7	Std.#8	Std.#9	Measurement Uncertainty (±°C)
35.0	34.76	34.86	34.81	34.83	35.22	35.21	35.22	35.12	34.83		0.28

Table 2: Reporting of Characterization Results

UUC* Setting (°C)	Max	UUC* Reading (°C)	Min	Average	Stability $\pm$ (°C)	Uniformity (°C)	Overall Variation (°C)
35.0	35.0	35.0	35.0	35.0	0.034	0.39	0.52

Note

UUC*: Unit Under Calibration

The quoted uncertainty includes contributions from Stability and Loading Effect (20% of Temperature Uniformity)

Stability = One-half of the greatest maximum difference of measured temperatures at any one sensor, for at least half an hour after reaching steady state

Uniformity = The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time

Overall Variation = The difference of The maximum and minimum measured temperatures throughout observation time

For more information, details of the sensor installation should be reported

From Air Chamber

Open Position

X Close 100% Fan-speed

Not Available

Dimension of inside chamber and sensor installation location (mm)

W x H x D: 65 x 125 x 48

W x H x D: 10 x 10 x 10



*The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k=2, providing a level of confidence of approximately 95%

End of Report

Signature

Calibration Certificate

Certificate No: 260233-006-01
Page: 1 of 2

Equipment: CHAMBER (Incubator)
Manufacturer: BINDER
Model: KB 400
Serial No: 2020000015530
ID No: UAE MIC 010/2564
Order No: 2602303
Operation No: 2602303 006
Date of Receipt: 11 March 2026
Date of Calibration: 11 March 2026

Customer Name: UNITED ANALYST AND ENGINEERING CONSULTANT CO., LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road, Bangkok, Phra Khanong, Bangkok 10260

Metrological traceability:

1. Calibration Method:

This instrument was calibrated by inserting 9 standard thermometers into its chamber, in accordance with W.T.E. 014 based on ILAS G.20.102.08 (E) Guidelines for Calibration and Checks of Temperature Controlled Enclosures. The temperature scale used in this laboratory is the International Temperature Scale of 1990 (ITS-90). All data shown below are final values and the final data may be obtained upon request.

2. Reference Standards Instrument

Instrument	Model/Type	Serial No./ID No	Certificate No.	Due Date	Through
Digital Thermometer with sensor	DAC61/3A	MY59004864	2502244-001-01	29 March 2026	NATIONAL FOOD INSTITUTE
	RTD	C-420-209 RTD0201-209			

3. This certificate is traceable to the International System of Units (SI)

4. This certificate is valid only for the instrument that was calibrated

5. The result of this calibration was found to be accurate only on the date and at the place of calibration shown

Calibrated by: Mr Nuttapol Niyomchat Specialist
Date of Issue: 18 March 2026

Approved by: (Mr Pheraphat Tuangit)
Manager, Division of Calibration Laboratory
Responsible for the Technical Management Team

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F: CS-009 Revision: 02 Date: 01/10/68

ศูนย์บริการและพัฒนาคุณภาพมาตรฐานอาหาร
ศูนย์บริการและพัฒนาคุณภาพมาตรฐานอาหาร
Foundation for Industrial Development National Food Institute
Food Industrial Laboratory Service Center

Certificate of Calibration

Certificate No: 260311-2-BL001-26
Code No: BL001-26
Page: 1 of 4

Customer Name: United Analyst and Engineering Consultant Co., Ltd.
Address: 3 Soi Udomsuk 41, Sukhumvit Rd., Bang Chak, Phra Khanong, Bangkok 10260

Equipment: Electronic Balance
Manufacturer: Mettler Toledo
Model: XSR205DU
Serial No: C210685394
Asset No: UAE WAO 010/2565

Building: UAE number 3 Floor: 2 Room: 208

Received Date: March 11, 2026

Date of Calibration: March 12, 2026

Calibration Conditions: Temperature: 23.0 °C to 23.6 °C
Humidity: 51.5 % to 54.4 %
Pressure: 755.7 mmHg to 757.9 mmHg

Calibrated by: Thanachat Chatpahol

Approved by: Suwit Chotnok

Signature: S. Chotnok

Issued Date: March 16, 2026

Note: 1) The Uncertainties are for a confidence probability of approximately 95%

2) This Certificate is valid only to the item calibrated on date and place of calibration

3) This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the United Analyst and Engineering Consultant Co., Ltd. (UAE)

Certificate No.: 260311-2-BL001-26
Code No.: BL001-26
Page: 2 of 4

Equipment: Electronic Balance
Model: XS205DU
Serial No.: C210685394
Max. Capacity: 81 g / 220 g
Calibration Date: March 12, 2026
Condition As-Received: in Condition

Manufacturer: Mettler Toledo
Readability: 0.1 to 81 (g) 0.01 mg, > 80 to 200 (g) 0.1 mg
ID No.: UAE WAO 010/2565

Condition of Equipment:
Condition of This Result of Calibration:
1. Calibration Method: This instrument was calibrated by method uM1 (F.6), 398 in-house Method based on JAS J40:2022
2. Reference Standards:
Reference Standard: Model Serial No. Calibrated By Certificate No. Traceability Due Date
Standard Weight Class 12 (DM1) 1 mg to 1 kg C310641605 METTLER TOLEDO GmbH 25-099359 METTLER TOLEDO GmbH 31 Mar 27
Standard Weight Class 12 (DM1) 1 mg to 1 kg B74919122 AMARC 25-099359 AMARC 21 Jan 27
Instrument Model Serial No. Calibrated By Certificate No. Traceability Due Date
Thermo Hygro Bare Meter RH-B 38253 Av. 86657 SUCCESS SCS-H 00992/68 Success Calibration 26 Nov 26
Thermo Hygro Bare Meter RH-B 38253 Av. 86657 TPA 18 Dec 26
3. This certification is traceable to SI unit.
4. This certification was certified only for the instrument we calibrated.
5. This result of calibration was found accurate as shown on date and place of calibration only.
6. Through the reference standard laboratory of Mettler-Toledo GmbH C310641605, SAS accreditation number SCS 0032
7. Through the reference standard laboratory of AMARC 25-099359, TSI accreditation number 0152
Calibration Result:
1. Repeatability of Reading:

Nominal Value (g)	Standard Deviation of Reading (g)
80	0.000019
200*	0.000055

2. Eccentric or off-center loading
A mass of 100 g was placed and moved to various position on pan.
The balance reading obtained is given in the table:

1	2	3	4	5	Maximum Difference (g)
(g)	(g)	(g)	(g)	(g)	
100.0002	100.0002	100.0003	100.0004	100.0003	0.0001

e.c.f
เอกสารไม่ควบคุม

Certificate No.: 260311-2-BL001-26
Code No.: BL001-26
Page: 3 of 4

Equipment: Electronic Balance
Model: XS205DU
Serial No.: C210685394
Max. Capacity: 81 g / 220 g
Calibration Date: March 12, 2026

Manufacturer: Mettler Toledo
Readability: 0.1 to 81 (g) 0.01 mg, > 80 to 200 (g) 0.1 mg
ID No.: UAE WAO 010/2565

Calibration Result: (Continued)
Calibration Range: 0 to 80 (g)
Calibration Adjustment: Internal Calibration

3. Error of indication from nominal or conventional mass value (Range: 0 to 81 (g); Resolution: 0.00001 g)					
Nominal Value (g)	Reference Value (g)	Indication (g)	Correction (g)	Uncertainty (± mg)	Coverage Factor k
Unload	0.0000000	0.00000	0.00000	0.032	2.65
0.001	0.0010020	0.00101	0.00000	0.032	2.65
0.005	0.0050025	0.00501	-0.00001	0.032	2.65
0.01	0.0100005	0.01001	-0.00001	0.032	2.65
0.05	0.0500003	0.05001	-0.00001	0.032	2.65
0.10	0.1000003	0.10000	0.00001	0.033	2.65
0.50	0.5000121	0.50001	0.00000	0.033	2.65
1	1.0000119	1.00003	-0.00002	0.032	2.52
2	2.0000125	2.00003	-0.00002	0.033	2.52
5	5.0000146	5.00003	-0.00002	0.033	2.52
10	10.0000330	10.00006	-0.00003	0.034	2.21
20	20.000028	20.00006	-0.00003	0.041	2.07
30	30.000061	30.00017	-0.00011	0.055	2.00
40	40.000054	40.00021	-0.00016	0.066	2.00
50	50.000037	50.00023	-0.00019	0.068	2.00
60	60.000070	60.00029	-0.00022	0.084	2.00
80	80.000098	80.00039	-0.00030	0.12	2.00

e.c.f

เอกสารไม่ควบคุม

Certificate No.: 260311-2-BL001-26
Code No.: BL001-26
Page: 4 of 4

Equipment: Electronic Balance
Model: XS205DU
Serial No.: C210685394
Max. Capacity: 81 g / 220 g
Calibration Date: March 12, 2026

Manufacturer: Mettler Toledo
Readability: 0.1 to 81 (g) 0.01 mg, > 80 to 200 (g) 0.1 mg
ID No.: UAE WAO 010/2565

Calibration Result: (Continued)
Calibration Range: > 81 to 200 (g)
Calibration Adjustment: Internal Calibration

3. Error of indication from nominal or conventional mass value (Range: > 81-220 (g); Resolution: 0.0001 g)					
Nominal Value (g)	Reference Value (g)	Indication (g)	Correction (g)	Uncertainty (± mg)	Coverage Factor k
Unload	0.0000000	0.0000	0.0000	0.032	2.65
90	90.0000910	90.0004	-0.0003	0.16	2.00
100	100.0000420	100.0003	-0.0003	0.16	2.00
110	110.0000750	110.0004	-0.0003	0.18	2.00
120	120.0000680	120.0004	-0.0003	0.19	2.00
130	130.0001010	130.0005	-0.0004	0.20	2.00
140	140.0000960	140.0004	-0.0003	0.21	2.00
150	150.0000790	150.0005	-0.0004	0.22	2.00
160	160.0001120	160.0004	-0.0003	0.23	2.00
170	170.0001050	170.0003	-0.0002	0.24	2.00
200	200.0001460	200.0001	0.0000	0.26	2.00

4. Effect of Tare test			
Tare Load (g)	Test Load (g)	Indication (g)	Correction (g)
100	20.000028	20.0001	0.0000
	40.000054	40.0001	0.0000
	60.000070	60.0000	0.0001
	80.000098	80.0000	0.0001
	100.000109	100.0000	0.0001

Remark: The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k, providing a level of confidence of approximately 95%.

0---0-End-0---0

e.c.f เอกสารไม่ควบคุม

Certificate of Calibration

Certificate No.: 260318-1-BL007-26
Code No.: BL007-26
Page: 1 of 3

Customer Name: United Analyst and Engineering Consultant Co., Ltd.
Address: 3 Soi Udomsuk 41, Sukhumvit Rd., Bang Chak, Phra Khanong, Bangkok 10260

Equipment: Electronic Balance
Manufacturer: Mettler Toledo
Model: AB204 S/FACT
Serial No.: 1129361010
Asset No.: UAE WAS 002/2552
Building: UAE No 3 Floor: 1 Room: 107
Received Date: April 2, 2026
Date of Calibration: April 3, 2026
Calibration Conditions: Temperature: start 21.9 °C end 22.1 °C
Humidity: start 60.6 % end 60.1 %
Pressure: start 756.6 mmHg end 755.1 mmHg

Calibrated by: Thanachat Chatpachol

Approved by: Suwit Chotnok

Signature: e.c.f

Issued Date: April 7, 2026

Note: 1) The Uncertainties are for a confidence probability of approximately 95%
2) This Certificate is valid only to the item calibrated on date and place of calibration.
3) This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme, which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the United Analyst and Engineering Consultant Co., Ltd. (LIAE)

เอกสารไม่ควบคุม

Equipment: Electronic Balance
 Model: AB204 S/FAC1
 Serial No.: 1129361010
 Max. Capacity: 220 g
 Calibration Date: Apr. 3, 2026
 Condition As-Received: in Condition

Manufacturer: Mettler Toledo
 Readability: 0.0001 g
 ID No.: UAE WAS 002/2552

Condition of Equipment

Condition of This Result of Calibration:

1. Calibration Method: This instrument was calibrated by the method of UAI (UAE) 004, 006.
2. Reference Standards: This instrument was calibrated by the method of UAI (UAE) 004, 006.
3. This certification is traceable to SI unit.
4. This certification was certified only for the instrument we calibrated.
5. This result of calibration was found accurate as shown on date and place of calibration only.
6. Through the reference standard laboratory of Mettler Toledo GmbH (S10641605, SAS accreditation number: SCS 0032).

Calibration Result:

1. Repeatability of Reading:

Nominal Value (g)	Standard Deviation of Reading (g)
200*	0.00010

2. Eccentric or off-center loading:

A mass of 100 g was placed and moved to various position on pan.
 The balance reading obtained is given in the table:

1	2	3	4	5	Maximum Difference (g)
(g)	(g)	(g)	(g)	(g)	
100.0025	100.0022	100.0022	100.0029	100.0030	0.0005

ei

เอกสารไม่ควบคุม

Equipment: Electronic Balance
 Model: AB204 S/FAC1
 Serial No.: 1129361010
 Max. Capacity: 220 g
 Calibration Date: Apr. 3, 2026

Manufacturer: Mettler Toledo
 Readability: 0.0001 g
 ID No.: UAE WAS 002/2552

Calibration Result: (Continued)
 Calibration Range: 0 - 200 g
 Calibration Adjustment: Internal Calibration

3. Error of indication from nominal or conventional mass value:

Nominal Value (g)	Reference Value (g)	Indication (g)	Correction (g)	Uncertainty ($\pm$ mg)	Coverage Factor k
Unloaded	0.0000000	0.0000	0.0000	0.17	2.37
0.01	0.0100005	0.0099	0.0001	0.17	2.32
0.05	0.0500003	0.0499	0.0001	0.17	2.32
0.1	0.1000083	0.0999	0.0001	0.17	2.32
0.5	0.5000121	0.5000	0.0000	0.17	2.32
1.0	1.0000119	1.0000	0.0000	0.17	2.32
10.0	10.000033	10.0003	-0.0003	0.17	2.32
40.0	40.000054	40.0010	-0.0010	0.18	2.23
50.0	50.000037	50.0013	-0.0012	0.18	2.23
80.0	80.000098	80.0020	-0.0019	0.19	2.11
100.0	100.000042	100.0025	-0.0024	0.20	2.10
120.0	120.000068	120.0028	-0.0027	0.22	2.06
150.0	150.000079	150.0035	-0.0035	0.24	2.04
160.0	160.000112	160.0036	-0.0035	0.25	2.03
200.0	200.000146	200.0049	-0.0047	0.28	2.00

4. Effect of Tare test:

Tare Load (g)	Test Load (g)	Indication (g)	Correction (g)
100	20.000028	20.0004	-0.0003
	40.000054	40.0009	-0.0008
	60.000070	60.0013	-0.0012
	80.000098	80.0017	-0.0016
	100.000042	100.0021	-0.0021

Remark:

The report uncertainty of measurement was found on standard uncertainty, multiplied by coverage factor k, providing a level of confidence of approximately 95%.

o-o-End-o-o-o

เอกสารไม่ควบคุม

List of Instrument Certificates for Environmental Quality Analysis

No.	Instrument/Equipment	Parameter	Manufacturer	Model/Serial No.	Calibrator	Certification No.	Date of Calibration	Due date of Calibration*
1	Analytical Balance	FAT OIL AND GREASE	Mettler Toledo	AB204-S/FACT / 1129361010	United Analyst and Engineering Consultant Co., Ltd.	250422 1 BL002 25	23/4/2025	22/4/2026
2	Analytical Balance	TOTAL DISSOLVED SOLIDS	Mettler Toledo	XSR205DU / C210685394	National Food Institute,Ministry of Industry, Thailand	25022226-002-01	20/3/2025	19/3/2026
3	Analytical Balance	TOTAL SUSPENDED SOLIDS	Mettler Toledo	XSR205DU / C009071872	National Food Institute,Ministry of Industry, Thailand	25022226-001-01	20/3/2025	19/3/2026
4	Auto Clave	TOTAL COLIFORM BACTERIA	ALP Co.,Ltd. (Japan)	CL-40L / 810010	National Food Institute Ministry of Industry (Thailand)	2503287-001-01	5/6/2025	4/6/2026
5	BOD Incubator	BIOCHEMICAL OXYGEN DEMAND	ARCO	UC4-1320 / 1021	Technology Promotion Association (Thailand-Japan)	24TM1002	7/7/2025	6/7/2026
6	DO Meter	BIOCHEMICAL OXYGEN DEMAND	YSI	5100 / 11B 101863	Technology Promotion Association (Thailand-Japan)	25TW29	17/2/2025	16/2/2026
7	Hot Air Oven	TOTAL DISSOLVED SOLIDS TOTAL SUSPENDED SOLIDS	Memmert	UF55 / B212.0411	Technology Promotion Association (Thailand-Japan)	25TM579	19/3/2025	18/3/2026
8	Cooled Incubator	TOTAL COLIFORM BACTERIA	Binder	KB400 / WTB20200000015535	National Food Institute, Ministry of Industry, Thailand	25022229-006-01	19/3/2025	18/3/2026
9	pH Meter	pH	Horiba	LAQUA-PH210 / HAOC0025	technology promotion association (thailand-japan	25CH261	26/2/2025	25/3/2026

Due Date of Calibration* : Based on the annual calibration plan. At least 1 time per year.



Certificate of Testing

Cert.No.: 25TW29
Page.: 1 of 2

Equipment : DO Meter
Manufacturer : YSI
Model : 5100
Serial No. : 11B 101863
ID No. : UAE.WAO.004/2554
Received Date : 14 February 2025
Test Date : 17 February 2025
Reference : 2502-0473DSC-1
Submitted by : United Analyst and Engineering Consultant Co., Ltd.
3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260
Laboratory Condition : Temperature (25 ± 5) °C
Humidity (50 ± 20) %
Test Procedure : In - house method : CP-CH9
by Comparison Technique with Azide Modification Method
Tested by : Walalak Sinithean
Approved by :
Approved Signatory
() Chakrit Waewwanjua
() Ponpan Paipim
(✓) Salthip Meangmai
Issue Date : 18 February 2025

เอกสารไม่ควบคุม



Cert.No.: 25TW29
Page.: 2 of 2

Condition of this result of calibration

1. Reference Standard Instruments :
This certification is traceable to the International System of Unit through the reference standards laboratory of Industrial Calibration Center, Technology Promotion Association (Thailand-Japan).

Instruments	Serial No.	ID No.	Certificate No.	Due Date
1. Burette	-	130BU10	23CG1172	22 Mar 2025
2. Balance	14233821	110RC001	24MM131	04 July 2025

2. Standard Material :-

Material	Manufacturer	Lot.No.	Assay
Sodium Thiosulfate 5-Hydrate AR	KEMAUS	2203162447	99.6%

Result : Dissolved Oxygen Meter Adjustment With Air 100 %
Dissolved Oxygen Probe No.: 24F100202

Titration Method (Azide Modification Method) (mg/L)	DO Meter Reading (mg/L)	Standard Deviation (mg/L)
8.22	8.22	0.0055

This report was certified only for the Instrument we tested it is allowable to use for study
Intend to use for advertising and referral purpose is prohibited. This report may not be reproduced
other in full, without written approval of the laboratory

-o-o-

เอกสารไม่ควบคุม



United Analyst and Engineering Consultant Co., Ltd.
3 Soi Udomsuk 41, Sukhumvit Road, Bangchak, Phrakhanong, Bangkok 10260
Tel: 02-2763-2828 Fax: 02-2763-2800 www.uaec consultant.com E-mail: uaec@uaec consultant.com



Certificate of Calibration

Certificate No.: 250422-1-BL002-25
Code No.: BL002-25
Page: 1 of 3

Customer Name: United Analyst and Engineering Consultant Co., Ltd.
Address: 3 Soi Udom suk 41, Sukhumvit Rd., Bang Chak, Phrakhanong, Bangkok 10260

Equipment: Electronic Balance
Manufacturer: Mettler Toledo
Model: AB204-S/FACT
Serial No.: 1129361010
Asset No.: UAE WAS 002/2552
Building : N/A Floor : 1 Room : 107

Received Date: April 22, 2025
Date of Calibration : April 23, 2025
Calibration Conditions: Temperature 22.8 °C to 23.4 °C
Humidity 54.8 % to 68.9 %
Pressure 756.6 mmHg to 758.2 mmHg

Calibrated by: Sakkarin Srirahang

Approved by: Suwit Chotnok

Signature:

Issued Date: April 25, 2025

- Note : 1) The Uncertainties are for a confidence probability of approximately 95%
2) This Certificate is valid only to the item calibrated on date and place of calibration
3) This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation Scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the United Analyst and Engineering Consultant Co., Ltd. (UAE)

เอกสารไม่ควบคุม

เอกสารไม่ควบคุม

Certificate No.: 250422-1-BL002-25
Code No.: BL002-25
Page: 2 of 3

Equipment: Electronic Balance
Model: AB204 S/FACT
Serial No.: 1129361010
Max. Capacity: 220 g
Calibration Date: April 23, 2025
Condition As-Received: In Condition

Manufacturer: Mettler Toledo
Readability: 0.0001 g
ID No.: UAE.WAS.002/2552

Condition of Equipment:

Condition of This Result of Calibration:

1. Calibration Method: This instrument was calibrated by method UAE.CP.CAL.006 In-House Method based on UNAS Lab 14 2022

2. Reference Standards:

Reference Standard:	Model	Serial No.	Calibrated By	Certificate No.	Traceability	Due Date
Standard Weight Class E2 (OIML)	1 mg to 1 kg	B749109122	AMARC	25 009359	Mettler Toledo	21 Jan 27
Standard Weight Class F1 (OIML)	1 mg to 200 g	11119512	AMARC	24 013840	Mettler Toledo	04 Feb 26
Instrument	Model	Serial No.	Calibrated By	Certificate No.	Traceability	Due Date
Thermo-Hydro-Baro Meter	MHB 3825D	AK-46457	SUCCESS	SG.H.00997/67	Success Gateway	21 Nov 25
Thermo-Hydro-Baro Meter	MHB 3825D	AK-46457	TPA	25P795	TPA	25 Feb 26

3. This certification is traceable to SI Unit

4. This certification was certified only for the instrument we calibrated

5. This result of calibration was found accurate as shown on date and place of calibration only

6. Through the reference standard laboratory of AMARC 25-009359 Calibration 0152

Calibration Result:


1. Repeatability of Reading:

Nominal Value (g)	Standard Deviation of Reading (g)
200*	0.00045

2. Eccentric or off-center loading

A mass of 100 g was placed and moved to various position on pan

The Balance reading obtained is given in the table

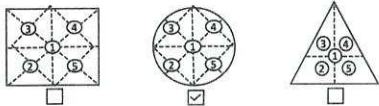

1 (g)	2 (g)	3 (g)	4 (g)	5 (g)	Maximum Difference (g)
100.0000	99.9996	99.9997	100.0003	100.0005	0.0005

เอกสารไม่ควบคุม

Calibration Certificate

Certificate No.: 2502229-006-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road,
Bangchack, Phrahanong, Bangkok 10260

Page 1 of 3

Equipment: CHAMBER (Incubator)

Manufacturer: BINDER

Model: KB 400

Serial No.: 20200000015535

ID No.: UAE.MIC.018/2564

Order No.: 2502229

Operation No.: 2502229-006

Date of Receipt: 19 March 2025

Date of Calibration: 19 March 2025

Calibrated by Mr.Jerawut Prapawuttipong
Scientist

Approved by *P. Jeyachandran*
(Mr.Pheraphat Tuanjit) (fr)

Manager, Division of Calibration Laboratory

Date of Issue: 25 March 2025 Responsible for the Technical Management Team

The uncertainties are for a confidence probability of approximately 95 %:

This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute

เอกสารไม่ควบคุม



Certificate No.: 250422-1-BL002-25
Code No.: BL002-25
Page: 3 of 3

Equipment: Electronic Balance
Model: AB204 S/FACT
Serial No.: 1129361010
Max. Capacity: 220 g
Calibration Date: April 23, 2025

Manufacturer: Mettler Toledo
Readability: 0.0001 g
ID No.: UAE.WAS.002/2552

Calibration Result: (Continued)

Calibration Range: 0 - 200 g

Calibration Adjustment: Internal Calibration

3. Error of indication from nominal or conventional mass value:

Nominal Value (g)	Reference Value (g)	Indication (g)	Correction (g)	Uncertainty (± mg)	Coverage Factor k
Unload	0.0000000	0.0000	0.0000	0.10	2.05
0.01	0.0100025	0.0099	0.0001	0.10	2.05
0.05	0.0500056	0.0500	0.0000	0.10	2.05
0.1	0.1000012	0.0999	0.0001	0.10	2.05
0.5	0.5000133	0.5000	0.0000	0.10	2.05
1	1.0000105	1.0000	0.0000	0.10	2.05
10	10.000010	10.0000	0.0000	0.11	2.04
40	40.000076	40.0000	0.0000	0.14	2.00
50	50.000056	50.0000	0.0001	0.13	2.00
80	80.000107	80.0000	0.0001	0.18	2.00
100	100.000109	99.9999	0.0002	0.17	2.00
120	120.00015	119.9999	0.0003	0.21	2.00
150	150.000165	149.9998	0.0003	0.24	2.00
160	160.000175	159.9997	0.0005	0.26	2.00
200	200.000129	199.9998	0.0004	0.30	2.00

4. Effect of Tare test:

Tare Load (g)	Test Load (g)	Indication (g)	Correction (g)
100	20.000041	19.9999	0.0001
	40.000076	39.9998	0.0002
	60.000066	59.9997	0.0003
	80.000107	79.9999	0.0002
	100.000168	100.0004	-0.0003

Remark:

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k, providing a level of confidence of approximately 95 %

เอกสารไม่ควบคุม

o--o-End-o--o

Calibration Report

Certificate No.: 2502229-006-01
Equipment: CHAMBER (Incubator)
Model: KB 400 Serial No.: 20200000015535
Resolution: 0.1 °C ID No.: UAE.MIC.018/2564
Manufacturer: BINDER

Date of Calibration: 19 March 2025 Page 2 of 3

Location: LABORATORY, UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.

Environment Condition: Ambient Temperature (18 ± 1) °C
Relative Humidity (50 ± 5) %
Line Voltage (223 ± 3) Volt

Condition of this results of Calibration:

1. This instrument was calibrated by insert 13 standard thermometer into its chamber and calibration according to W-TE-014 Based on TLAS G 20-1/02-08 (E) Guidelines for Calibration and Checks of Temperature Controlled Enclosures.
The temperature scale used was based on ITS - 90
All data show below were final values and the initial data may be obtained upon request.

2. Reference Standard Instrument:

Instrument	Model	Serial No./ID No.	Certificate No.	Due Date	Through
Digital Thermometer with sensor	34972A	MY49616851	TE 670477-01	4 May 2025	NATIONAL FOOD INSTITUTE
	RTD	CH-4961303 / H10-201303			

3. This certificate is traceable to International System of Units (SI Units)

4. This certificate was certified only for the instrument we calibrated

5. This result of calibration was found accurate as shown on date and place of calibration only

6. Condition of Calibrated item: Good

UUC Description:

Time of Record: 1 Hour 9 Minute At: 35.0 °C
Fresh air Damper: - Open Position: -
X Close Fan: -
Not Available

7. Result of Calibration: X Without adjustment After adjustment

เอกสารไม่ควบคุม



Calibration Report

Certificate No.: 2502229-006-01
Equipment: CHAMBER (Incubator)
Model: KB 400 Serial No.: 2020000015535
Resolution: 0.1 °C ID No.: UAE MIC 018/2564
Manufacturer: BINDER

Date of Calibration: 19 March 2025 Page 3 of 3

Calibration point: 35.0 °C

Calibration result:

Calibration Condition	Temperature (°C)	Relative Humidity (%)	Line Voltage (VOLT)
MIN	17.1	45	220.0
MAX	18.1	55	225.0

Table 1 : Reporting of Temperature

Calibration point (°C)	Measured Temperature (°C) @ Sensor No. (Sensor No.13 is REF)													Uncertainty ± (°C)
	#1	#2	#3	#4	#5	#6	#7	#8	#9	#10	#11	#12	#13	
35.0	34.98	35.17	34.99	34.92	35.18	35.01	35.00	35.13	35.00	34.96	35.02	35.17	35.04	0.27

Table 2 : Reporting of Characterization Result

UUC* Setting (°C)	UUC* Reading (°C)			Temperature Stability ± (°C)	Temperature Uniformity (°C)	Overall Variation (°C)
	MIN	MAX	Average			
35.0	35.0	35.0	35.0	0.029	0.15	0.30

Note: The quoted uncertainty include "Stability" and "Loading effect (20% of Temp Uniformity)"

UUC* = Unit Under Calibration

Stability = One-half of the greatest maximum difference of measured temperatures at any one sensors, for at least half an hour after reaching steady state

Uniformity = The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time

Overall Variation = The difference of the maximum and minimum measured temperatures throughout observation time

The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k= 2, providing a level of confidence of approximately 95%

----- End -----

F: CS-012 Revision: 01 Date: 20-04-65

เอกสารไม่ควบคุม
25 March 2025

2503 0437 OC-3
2503 0437 OC-3
2503 0437 OC-3
2503 0437 OC-3

Certificate of Calibration

Cert. No.: 25TM579

Page : 1 of 3

Equipment : Hot Air Oven

Manufacturer : Memmert

Model : UF 55

Serial No. : B212.0411

ID No. : UAE WAO.005/2556

Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260

Location : Lab Floor 2

Received Order : 19 March 2025

Calibration Date : 19 March 2025

Ambient Temperature : (26 ± 10) °C

Relative Humidity : (50 ± 30) %

AC Line Voltage : (220 ± 22) V

Calibrated by : Man Pattanapongpaiboon

Approved by :

Kunchit

Approved Signatory

() Chakrit Waewwanjua

() Suwit Imjai

(✓) Kunchit Promprat

Issue Date : 27 March 2025

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written approval of the head of Corporate Services 3: Equipment Calibration and Testing Services

เอกสารไม่ควบคุม



Equipment : Hot Air Oven
Condition As-Received : Used Item
Reference : 2503-0437OC-3

Cert. No.: 25TM579
Page : 2 of 3

Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD) and Thermocouple Type T.

The temperature scale used was based on ITS-90.

Condition of this result of calibration

- Reference standard instrument-
Instrument Serial No. Cert. No. Traceable Due Date
1) Data Acquisition MY44073381 24LM73 TPA 18 May 2025
- This certificate is valid only to the item calibrated on date and place of calibration.
- This certificate is traceable to the International System of Unit.

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration :- (*) Without Adjustment


Function of UUC* : Temperature Source

Fresh air setting : Close

Environment during calibration	
	Beginning
Temp. (°C)	27
REL Humid. (%)	49
AC Supply (Volt)	221
	224

Ref. Std. ID No.: @ Calibration Point

Position :	(120,180) °C	(104) °C
1	23-01TC-01	1RTD-2/1
2	23-01TC-02	1RTD-2/2
3	23-01TC-03	22-01RTD-03
4	23-01TC-04	1RTD-2/4
5	23-01TC-05	1RTD-2/5
6	23-01TC-06	1RTD-2/6
7	23-01TC-07	23-01RTD-07
8	23-01TC-08	1RTD-2/8
9 (ref.)	23-01TC-09	23-01RTD-09

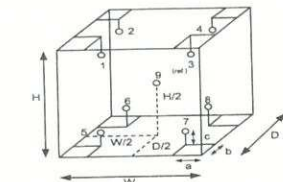

Probe Installation Details : Dimension of Chamber :
a = 5.0 cm D = 0.50 m
b = 5.0 cm W = 0.80 m
c = 5.0 cm H = 0.75 m
Capacity = 0.30 m³



Equipment : Hot Air Oven
Condition As-Received : Used Item
Reference : 2503-0437OC-3

Cert. No.: 25TM579
Page : 3 of 3

Result of Calibration :-

Function of UUC* : Temperature Source

Fresh air setting : Close

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor k
104.0	104.0	104.0	0.040	0.43	0.78	2
120.0	120.0	120.0	0.64	1.3	1.6	2
180.0	180.0	180.0	0.49	1.5	1.8	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (±°C)
	1	2	3	4	5	6	7	8	9 (ref.)	
104.0	104.335	104.135	104.363	104.317	103.649	103.738	104.179	104.229	104.025	0.42
120.0	119.575	119.366	119.807	119.905	118.994	119.194	119.888	119.994	120.064	1.1
180.0	180.286	179.510	180.401	180.551	179.281	179.463	180.196	180.451	180.374	1.2

Average* : The average of 30 values in each position.

Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor

Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.

Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation.

UUC* : Unit Under Calibration

Note : The reported uncertainty of measurement was included stability and excluded uniformity .

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่ควบคุม

เอกสารไม่ควบคุม

Calibration Certificate

Certificate No.: 2503287-001-01
Client name: UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Address: 3 Soi Udomsuk 41, Sukhumvit Road,
Bangchack, Prakhonong, Bangkok 10260

Page 1 of 3

Equipment: Autoclave
Manufacturer: ALP
Model: CL-40L
Serial No.: 810010
ID No.: UAE.MIC.032/2565
Order No.: 2503287
Operation No.: 2503287-001
Date of Receipt: 5 June 2025
Date of Calibration: 5 June 2025

Calibrated by Mr.Pheraphat Tuanjit
Scientist
Approved by *P. Jenghabit*
(Miss Preeyaporn Jaengkarnkit)
Vice President, Department of Laboratory Services
Responsible for the Technical Management Team
Date of Issue: 11 June 2025

The uncertainties are for a confidence probability of approximately 95 %.
This Certificate is issued in accordance with the conditions of accreditation granted by the Thai Laboratory Accreditation scheme which has assessed the measurement capability of the laboratory and its traceability to recognized national standards and to the units of measurement realized at the corresponding national standards laboratory. This certificate may not be reproduced other than in full except with the prior written approval of the National Food Institute.

F-CS-009 Revision: 01 Date: 20-04-65

2503287-001-01 : Autoclave Calibration Certificate
2503287-001-01 : Autoclave Calibration Certificate
2503287-001-01 : Autoclave Calibration Certificate
Tel : +66(0) 2-422 6555 Fax : +66(0) 2-422 6555

Calibration Report

Certificate No.: 2503287-001-01
Equipment: Autoclave
Model: CL-40L
Resolution: 1 °C
Manufacturer: ALP
Serial No.: 810010
ID No.: UAE.MIC.032/2565

Date of Calibration: 5 June 2025
Calibration point: 115 and 121 °C

Calibration Condition	Temperature (°C)	Relative Humidity (%)	Line Voltage (Volt)
Min	25.8	50	225
Max	26.8	59	235

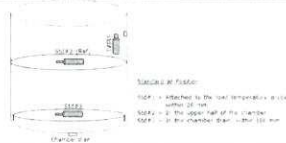


Table 1 : Reporting of Temperature

Calibration Point (°C)	Measured Temperature (°C) @ Sensor No. (Sensor No.2 is REF)			Uncertainty ± (°C)
	Std.# 1	Std.# 2 (Ref)	Std.# 3	
115	115.46	115.43	115.42	0.70
121	121.59	121.54	121.51	0.70

Table 2 : Reporting of Characterization Result

UUC* Setting (°C)	UUC* Reading			Stability ± (°C)	Uniformity (°C)	Overall Variation (°C)
	Min (°C)	Max (°C)	Average (°C)			
115	115	115	115	0.08	0.17	0.50
121	121	121	121	0.12	0.19	0.52

NOTE

The quoted uncertainty include " Stability " and " Loading effect (20% of Uniformity)"
UUC* = Unit Under Calibration
Stability = One-half of the greatest maximum difference of measured temperatures at any one sensors, for at least half an hour after reaching steady state.
Uniformity = The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time.
Overall Variation = The difference of the maximum and minimum measured temperatures throughout observation time.
The report uncertainty of measurement was based on standard uncertainty multiplied by coverage factor k= 2, providing a level of confidence of approximately 95 %.

----- End -----

F-CS-012 Revision: 01 Date: 20-04-65

2503287-001-01 : Autoclave Calibration Certificate
2503287-001-01 : Autoclave Calibration Certificate
2503287-001-01 : Autoclave Calibration Certificate
Tel : +66(0) 2-422 6555 Fax : +66(0) 2-422 6555

Calibration Report

Certificate No.: 2503287-001-01
Equipment: Autoclave
Model: CL-40L
Resolution: 1 °C
Manufacturer: ALP
Serial No.: 810010
ID No.: UAE.MIC.032/2565
Date of Calibration: 5 June 2025

Page 2 of 3

Location: Room 301 Media Preparation, UNITED ANALYST AND ENGINEERING CONSULTANT CO.,LTD.
Environment Condition: Ambient Temperature (26 ± 1) °C
Relative Humidity (55 ± 5) %
Line Voltage (230 ± 5) Volt

Condition of this results of Calibration:

- This instrument was calibrated by insert 3 standard Data loggers with RTD into its autoclave and calibration according to W-TE-018 based on BS 2646-1:2021, Autoclaves for sterilization in laboratories
Part 1: Design, construction, safety and performance - Specification
- The temperature scale used was based on ITS - 90.
- All data show below were final values and the initial data may be obtained upon request.
- Reference Standard Instrument :

Instrument	Model	Serial No.	Certificate No.	Due Date	Through
Digital Thermometer with RTD (Data Logger)	HiTemp140-PT	T20627	NC-25-03-18-181	11 Mar-26	MADETECH, INC.
	QM-CP-HITEMP-140	R56916	2502081-002-01	11 Mar-26	NATIONAL FOOD INSTITUTE
	PRTemp140	R38546	2501835-001-01	22 Feb-26	NATIONAL FOOD INSTITUTE

- This certificate is traceable to International System of Units (SI Units).
- This certificate was certified only for the instrument we calibrated.
- This result of calibration was found accurate as shown on date and place of calibration only.
- This standard does not apply to sterilizers or disinfectors used for medical, dental, pharmaceutical.
- Condition of Calibrated item : Good
UUC Description : Setting program function sterilization : STERILIZE/NORMAL
Time of sterilization : 20 Minute At : 115 and 121 °C
- Result of Calibration : ☒ Without adjustment
☐ After adjustment

P. Jenghabit
11 June 2025

F-CS-012 Revision: 01 Date: 20-04-65

2503287-001-01 : Autoclave Calibration Certificate
2503287-001-01 : Autoclave Calibration Certificate
2503287-001-01 : Autoclave Calibration Certificate
Tel : +66(0) 2-422 6555 Fax : +66(0) 2-422 6555



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL 0-2717-3000-29 FAX 0-2719-9484

Certificate of Calibration

Cert.No.: 25CH261
Page: 1 of 3

Equipment : pH Meter
Manufacturer : Horiba
Model : LAQUA-PH210
Serial No. : HA0C0025
ID No. : UAE.EFM.117/2563(EFM.pH.07/63)

Condition As-Received: Used Item
Received Date : 25 February 2025
Calibration Date : 26 to 28 February 2025
Reference : 2502-0783WSC-1
Submitted by : United Analyst and Engineering Consultant Co.,Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchack, Phrakhanong, Bangkok 10260

Ambient Temperature : (25 ± 2.5) °C
Relative Humidity : (50 ± 15) %
Calibration Procedure : In - house method :
- CP-CH5 by direct measurement with DC voltage standard and direct measurement with certified reference material (CRM)
- CP-CH6 by comparison with temperature standard

Calibrated by : Warakorn Lemagatrakul

Approved by : *P. Jenghabit*
Approved Signatory

() Chakrit Waewwanjua
() Ponpan Paipim
(✓) Saithip Meangmai

Issue Date : 28 February 2025

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written Approval of the head of Corporate Services 3: Equipment Calibration and Testing Services.

เอกสารไม่ควบคุม



Cert.No.: 25CH261
Page.: 2 of 3

Condition of this calibration result

1. Reference Standard Instrument

Instrument	Serial No.	ID No.	Cert. No.	Due Date
1) Document Process Calibrator	54030049	130RC116	24E2759	25 Aug 2025
2) Ref. Standard Thermometer	4982054	110RC044	24I757	14 July 2025

- This Certification is traceable to SI Through Technology Promotion Association (Thailand - Japan)

2. Certified Reference Materials : The measurement results are traceable to SI through Hach Lenge GmbH Ltd., Deutsche Akkreditierungsstelle, Accredited No. D-RM-15184-01-00
: The measurement results are traceable to SI through CPA chem Ltd., ANSI-ASQ National Accreditation Board, Accredited No. AR-1835

Buffer Solution	Manufacturer	Lot No.	Exp. date
pH 4.007	CPA chem	1066665	18 Jan 2027
pH 6.999	Hach Lenge GmbH	C03220	29 Oct 2026
pH 10.010	CPA chem	1066669	18 Jan 2026

3. This certificate is valid only to the item calibrated on date and place of calibration.

Calibration Results

Function : mV Measurement

Performing standard curve by Document Process Calibrator at pH (4,7)(7,10)

Unit Under Calibration	Nominal Value	Standard Voltage Input	Actual Reading		Uncertainty of Measurement ($\pm$ mV)	Coverage factor k
	pH	mV	mV	pH		
pH Meter S/N.: HAOC0025	4.00	177.48	177.5	4.01	0.058	2.00
	7.00	0.00	0.0	7.02	0.058	2.00
	7.00	0.00	0.0	7.02	0.058	2.00
	10.00	-177.48	-177.5	10.01	0.058	2.00



Cert.No.: 25CH261
Page.: 3 of 3

Calibration Results

Function : pH Measurement

Performing three buffers standard curve by using buffer nominal pH (4,7)(7,10)

Unit Under Calibration	Standard pH Buffer Solution	Actual pH Reading	Actual mV Reading (mV)	Uncertainty of pH Measurement ($\pm$)	Coverage factor k
pH Electrode S/N.: Q9AG0214	4.007	4.01	178.4	0.0071	2.00
	6.999	7.00	4.1	0.0092	2.00
	6.999	7.00	3.0	0.0095	2.00
	10.010	10.01	-169.8	0.0092	2.00

Function : Temperature Measurement

(*) Without adjustment

This equipment was connected with Temperature Probe;

- Model : 9652-10D
- Serial No. : Q9AG0214

Dimension of probe

- Length : 110 mm.
- Diameter : 16 mm.
- Immersion Depth : 80 mm.

Calibration Point (°C)	Standard Temperature (°C)	UUC* Reading (°C)	Error (°C)	Uncertainty of measurement ($\pm$ °C)	Coverage factor k
15.0	15.002	15.0	-0.002	0.13	2.00
30.0	30.003	30.0	-0.003	0.13	2.00
45.0	45.002	44.9	-0.102	0.13	2.00

Remark : - UUC* = Unit Under Calibration

The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor k, providing a level of confidence of approximately 95 %.

-o-o-

เอกสารไม่ควบคุม

เอกสารไม่ควบคุม



TECHNOLOGY PROMOTION ASSOCIATION (THAILAND-JAPAN)
CORPORATE SERVICES 3: EQUIPMENT CALIBRATION AND TESTING SERVICES
534/4 PATTANAKARN ROAD SOI 18, SUANLUANG, SUANLUANG BANGKOK 10250
TEL 0-2717-3000-29 FAX 0-2719-9484



Certificate of Calibration

Cert. No.: 25TM1002
Page : 1 of 3

Equipment : BOD Incubator
Manufacturer : ARCO
Model : UC4-1320
Serial No. : -
ID No. : UAE.WAO.002/2550
Submitted by : United Analyst and Engineering Consultant Co., Ltd.
3 Soi Udomsuk 41, Sukhumvit Road,
Bangchak, Phrakhanong,
Bangkok 10260
Location : Lab Floor 2
Received Order : 07 July 2025
Calibration Date : 07 July 2025
Ambient Temperature : (26 $\pm$ 10) °C
Relative Humidity : (50 $\pm$ 30) %
AC Line Voltage : (220 $\pm$ 22) V

Calibrated by : Man Pattanapongpaiboon

() Chakrit Waewwanjua
() Suwit Imjai
(✓) Kunchit Promprat

Issue Date : 17 July 2025

The Uncertainties are for a confidence probability of approximately 95%

This certificate may not be reproduced other than in full, except with the prior written
Approval of the head of Corporate Services 3: Equipment Calibration and Testing Services

เอกสารไม่ควบคุม



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2507-0146OC-3

Cert. No.: 25TM1002
Page : 2 of 3

Procedure Used :-

Calibration were conducted using calibration procedure CP-OT02 based on TLAS G-20 according to direct measurement method with Data Acquisition which connected with Resistance Temperature Detector (RTD).
The temperature scale used was based on ITS-90.

Condition of this result of calibration

1. Reference standard instrument:-

Instrument	Serial No.	Cert. No.	Traceable	Due Date
1) Data Acquisition	MY59003411	24LM192	TPA	24 Dec 2025

2. This certificate is valid only to the item calibrated on date and place of calibration.


3. This measurement result is traceable to the International System of Unit maintained through :

Remark : TPA : Technology Promotion Association (Thailand - Japan)

Result of Calibration : (*) Without Adjustment

Function of UUC* : Temperature Source

Fresh air setting : Not Available

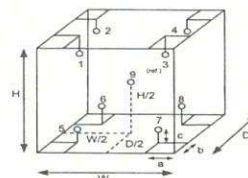

Probe Installation Details :

a = 10 cm
b = 10 cm
c = 10 cm

Dimension of Chamber :

D = 0.62 m
W = 1.2 m
H = 1.2 m
Capacity = 0.89 m³

Environment during calibration		
	Beginning	Finished
Temp. (°C)	28	28
REL.Humid. (%)	57	58
AC Supply (Volt)	225	226

Position :	Ref. Std. ID No.
1	25-20RTD-2/1
2	25-20RTD-2/2
3	25-20RTD-2/3
4	20RTD-2/4
5	20RTD-2/5
6	20RTD-2/6
7	20RTD-2/7
8	20RTD-2/8
9 (ref.)	20RTD-2/9

เอกสารไม่ควบคุม



Equipment : BOD Incubator
Condition As-Received : Used Item
Reference : 2507-0146OC-3
Result of Calibration :-
(*) Without Adjustment
Function of UUC* : Temperature Source
Fresh air setting : Not Available

Cert. No.: 25TM1002
Page : 3 of 3

Calibration Point (°C)	UUC* Setting (°C)	UUC* Reading (°C)	Temperature stability (± °C)	Temperature uniformity (°C)	Overall Variation (°C)	Coverage Factor <i>k</i>
20.0	20.0	19.8	0.48	0.46	1.2	2

Calibration Point (°C)	Measured Temperature (°C)									Uncertainty (± °C)
	Position									
	1	2	3	4	5	6	7	8	9 (ref.)	
20.0	20.226	20.266	20.121	19.807	19.971	20.055	19.872	19.853	19.925	0.72

Average* : The average of 30 values in each position.
Temperature stability : One-half of the greatest maximum difference of measured temperature at any one sensor.
Temperature uniformity : The maximum difference of measured temperatures at any sensors and the measured temperature at the reference location which are observed at the same time or at as close an observation time as possible to determine the temperature pattern or homogeneity within the chamber under steady-state conditions.
Overall Variation : The Difference of the maximum and minimum measured temperatures throughout observation
UUC* : Unit Under Calibration
Note : The reported uncertainty of measurement was included stability and excluded uniformity .
The reported uncertainty of measurement was based on a standard uncertainty multiplied by a coverage factor *k*, providing a level of confidence of approximately 95 %.

-o0o-

เอกสารไม่ควบคุม

ภาคผนวก จ
หนังสือรับรองการขึ้นทะเบียน
ห้องปฏิบัติการวิเคราะห์เอกชน

ที่ ออ ๐๓๐๘/๒๕๖๓



กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๒ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๕๐๐

๒๕ สิงหาคม ๒๕๖๔

เรื่อง เอกสารแนบท้ายหนังสือเปลี่ยนแปลงบุคลากรและสารมลพิษที่วิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท ยูโนเดค แอนนาลิติกส์ แอนด์ เอ็นจิเนียริ่ง คอนซัลแตนท์ จำกัด

ปิ่นแก้ว ศาลาพระที่นั่ง/ตึกอยู่/เปลี่ยนแปลงบุคลากร และขอแจ้งสารมลพิษที่ขอปฏิบัติวิเคราะห์แยกตาม
ลงวันที่ ๒ สิงหาคม ๒๕๖๔

ตามที่ขอที่แจ้ง บริษัท ยูโนเดค แอนนาลิติกส์ แอนด์ เอ็นจิเนียริ่ง คอนซัลแตนท์ จำกัด
ห้องปฏิบัติการวิเคราะห์แยกตาม เลขทะเบียน ๖-๑๕๕ สถานที่ตั้งเลขที่ ๑ ซอยอุดมสุข ๔๑ ถนนสุขุมวิท
แขวงบางจาก เขตพระโขนง กรุงเทพมหานคร ขอแจ้งบุคลากร และขอแจ้งสารมลพิษที่ขอปฏิบัติวิเคราะห์แยกตาม
ลงวันที่ ๒ สิงหาคม ๒๕๖๔

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว ให้ยกเลิกเจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์แยกตาม
จำนวน ๕ ราย ได้แก่

- | | |
|-----------------------------|----------------------------|
| ๑) น.ช.ณัฐดนัย พงษ์ศิริธรรม | ทะเบียนเลขที่ ๖-๑๕๕-๑-๐๑๐๔ |
| ๒) นายธีรวัฒน์ ธรรมสุวรรณ | ทะเบียนเลขที่ ๖-๑๕๕-๑-๐๑๐๕ |
| ๓) นายอาทิตย์ ส.ภา | ทะเบียนเลขที่ ๖-๑๕๕-๑-๐๑๐๖ |
| ๔) นางสาวกนกชนก ภูมิคำ | ทะเบียนเลขที่ ๖-๑๕๕-๑-๐๑๐๗ |
| ๕) นายวิรัชพงศ์ แสงทำนุ | ทะเบียนเลขที่ ๖-๑๕๕-๑-๐๑๐๘ |

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

(นายศิระ จันทร์โอ)

ผู้อำนวยการกองบริหารสุขภาพ
ผู้ชำนาญการพิเศษและผู้อำนวยการกอง
ปฏิบัติการทางอนามัยสิ่งแวดล้อม

กองวิจัยและเตือนภัยมลพิษทางอากาศ

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและประเมินภัยมลพิษทางอากาศ

โทร. ๐ ๒๕๓๐ ๒๕๓๖ ต่อ ๒๐๑๕๕

โทรสาร ๐ ๒๕๓๐ ๒๕๓๖ ต่อ ๒๐๑๕๕

ไปรษณีย์อิเล็กทรอนิกส์ sarabangdw@mail.go.th



ดำเนินการถูกต้อง



ผู้ชำนาญการพิเศษและผู้อำนวยการกองบริหารสุขภาพ

ที่ ออ ๐๓๐๘/๒๕๖๓



กรมโรงงานอุตสาหกรรม
ถนนพระรามที่ ๒ แขวงทุ่งพญาไท
เขตราชเทวี กรุงเทพฯ ๑๐๕๐๐

๐๗ กรกฎาคม ๒๕๖๔

เรื่อง เปลี่ยนแปลงบุคลากรและสารมลพิษที่วิเคราะห์

เรียน กรรมการผู้จัดการ บริษัท ยูโนเดค แอนนาลิติกส์ แอนด์ เอ็นจิเนียริ่ง คอนซัลแตนท์ จำกัด

อ้างถึง คำขอขึ้นทะเบียน/ตึกอยู่/เปลี่ยนแปลงบุคลากร และขอแจ้งสารมลพิษที่ขอปฏิบัติวิเคราะห์แยกตาม
ลงวันที่ ๒ สิงหาคม ๒๕๖๔

สิ่งที่ส่งมาด้วย เอกสารแนบท้ายหนังสือเปลี่ยนแปลงบุคลากรและสารมลพิษที่วิเคราะห์

บริษัท ยูโนเดค แอนนาลิติกส์ แอนด์ เอ็นจิเนียริ่ง คอนซัลแตนท์ จำกัด จำนวน ๒ ฉบับ

ตามที่ขอที่แจ้ง บริษัท ยูโนเดค แอนนาลิติกส์ แอนด์ เอ็นจิเนียริ่ง คอนซัลแตนท์ จำกัด
ห้องปฏิบัติการวิเคราะห์แยกตาม เลขทะเบียน ๖-๑๕๕ สถานที่ตั้งเลขที่ ๑ ซอยอุดมสุข ๔๑ ถนนสุขุมวิท
แขวงบางจาก เขตพระโขนง กรุงเทพมหานคร ขอแจ้งบุคลากร และขอแจ้งสารมลพิษที่ขอปฏิบัติวิเคราะห์
ความเปลี่ยนแปลงดังนี้

กรมโรงงานอุตสาหกรรมพิจารณาแล้ว มีความเห็นดังนี้

๑. โดยยกเลิกเจ้าหน้าที่ห้องปฏิบัติการวิเคราะห์แยกตาม จำนวน ๓ ราย

๑) นายสุเชษฐ์ พันสิงห์ ทะเบียนเลขที่ ๖-๑๕๕-๑-๐๑๐๑

๒) นางสาวสุกัญญา เอี่ยมเงิน ทะเบียนเลขที่ ๖-๑๕๕-๑-๐๑๐๒

๓) นางสาวชานันดา กีมาคม ทะเบียนเลขที่ ๖-๑๕๕-๑-๐๑๐๓

๒. ให้เพิ่มผู้ควบคุมห้องปฏิบัติการวิเคราะห์แยกตาม จำนวน ๓ ราย

นายสุเชษฐ์ พันสิงห์ ทะเบียนเลขที่ ๖-๑๕๕-๑-๐๑๐๑

๓. ให้เพิ่มสายส่งสารมลพิษที่วิเคราะห์ในน้ำได้ดิน ยากาเคียว และดิน ดินสิ่งส่งมาด้วย

อนึ่ง หนังสือฉบับนี้จะแนบมาพร้อมหนังสือรับทราบเปลี่ยนแปลงห้องปฏิบัติการวิเคราะห์แยกตาม

ในวันที่ ๒ กุมภาพันธ์ ๒๕๖๔

จึงเรียนมาเพื่อทราบ

ขอแสดงความนับถือ

(นายศิระ จันทร์โอ)

ผู้อำนวยการกองบริหารสุขภาพ
ผู้ชำนาญการพิเศษและผู้อำนวยการกอง
ปฏิบัติการทางอนามัยสิ่งแวดล้อม

กองวิจัยและเตือนภัยมลพิษทางอากาศ

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและประเมินภัยมลพิษทางอากาศ

โทร. ๐ ๒๕๓๐ ๒๕๓๖ ต่อ ๒๐๑๕๕

โทรสาร ๐ ๒๕๓๐ ๒๕๓๖ ต่อ ๒๐๑๕๕

ไปรษณีย์อิเล็กทรอนิกส์ sarabangdw@mail.go.th



ดำเนินการถูกต้อง



ผู้ชำนาญการพิเศษและผู้อำนวยการกองบริหารสุขภาพ

เอกสารแนบท้ายหนังสือเปลี่ยนแปลงบุคลากรและสารมลพิษที่วิเคราะห์

บริษัท ยูโนเดค แอนนาลิติกส์ แอนด์ เอ็นจิเนียริ่ง คอนซัลแตนท์ จำกัด เลขทะเบียน ๖-๑๕๕

ที่ ออ ๐๓๐๘/๒๕๖๓

ลงวันที่ ๐๗ กรกฎาคม ๒๕๖๔

ขอแจ้งสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๑๑ รายการ

น้ำได้ดิน จำนวน ๔ รายการ

ลำดับที่	สารมลพิษ	วิธีการวิเคราะห์
1	Aluminum	Digestion, Inductively Coupled Plasma Method ⁽¹⁾
2	Copper	1) Digestion, Flame Atomic Absorption Spectrometric Method ⁽¹⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽¹⁾
3	Iron	1) Digestion, Flame Atomic Absorption Spectrometric Method ⁽¹⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽¹⁾
4	Molybdenum	Digestion, Inductively Coupled Plasma Method ⁽¹⁾

อากาศเคียว (ปล่องระเหย) จำนวน ๑ รายการ

ลำดับที่	สารมลพิษ	วิธีการวิเคราะห์
1	Oxides of Nitrogen	Absorption Sampling, Ion Chromatographic Method ⁽¹⁾

ดิน จำนวน ๖ รายการ

ลำดับที่	สารมลพิษ	วิธีการวิเคราะห์
1	Aluminum	Digestion, Inductively Coupled Plasma Method ^(3,5)
2	Copper	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(3,6) 2) Digestion, Inductively Coupled Plasma Method ^(3,5)
3	Iron	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(3,6) 2) Digestion, Inductively Coupled Plasma Method ^(3,5)
4	Molybdenum	Digestion, Inductively Coupled Plasma Method ^(3,6)
5	pH	Electrometric Method ⁽⁸⁾
6	TPH (C ₂₅ -C ₆)	Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(4,7)

เอกสารอ้างอิง

1. APHA, AWWA, WEF. Standard Methods for the Examination of Water and Wastewater. 24th ed. Washington, DC: APHA, 2023.

2. United States Environmental Protection Agency. Standards of Performance for New Stationary Sources. 40 CFR 60. Appendix A, 2023.

3. United States...

กลุ่มมาตรฐานวิธีการวิเคราะห์ทดสอบมลพิษและประเมินภัยมลพิษทางอากาศ การวิจัยและเตือนภัยมลพิษทางอากาศ กรมโรงงานอุตสาหกรรม โทร. ๐ ๒๕๓๐ ๒๕๓๖ ต่อ ๒๐๑๕๕

เอกสารแนบท้ายหนังสือตอบข้อหารือเรื่องระเบียบห้องปฏิบัติการวิเคราะห์เอกชน
บริษัท ปูนซีเมนต์ แอนนาทิสต์ แอนด์ เอ็นจิเนียริง คอนสตรัคชั่น จำกัด เลขทะเบียน ๖-๔๔๙
ที่ ออก ๐๓๐๐(๑)๗ ๑๐ ๔ ๗ ลงวันที่ ๐๗ กุมภาพันธ์ ๒๕๖๕

ขอแจ้งสารมลพิษที่ได้รับขึ้นทะเบียนจากกรมโรงงานอุตสาหกรรม จำนวน ๓๘๗ รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Aldrin	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
2	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
3	Barium	Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
4	α-BHC	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
5	β-BHC	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
6	δ-BHC	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
7	γ-BHC	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
8	Biochemical Oxygen Demand	1) 5-Day BOD Test, Azide Modification Method ⁽⁴⁾ 2) 5-Day BOD Test, Membrane Electrode Method ⁽⁴⁾
9	Cadmium	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
10	Chemical Oxygen Demand	1) Closed Reflux, Titrimetric Method ⁽⁴⁾ 2) Closed Reflux, Colorimetric Method ⁽⁴⁾ 3) Open Reflux, Titrimetric Method ⁽⁴⁾
11	Chlordane	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
12	Chromium	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
13	Color	ADMI Weighted-Ordinate Spectrophotometric Method ⁽⁴⁾
14	Copper	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
15	Cyanide	1) Distillation, Colorimetric Method ⁽⁴⁾ 2) Total Cyanide after Distillation, by Flow Injection Analysis Method ⁽⁴⁾
16	o,p'-DDT	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
17	4,4'-DDD	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
18	4,4'-DDE	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
19	4,4'-DDT	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
20	Dieldrin	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
21	Endosulfan I	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
22	Endosulfan II	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
23	Endosulfan sulfate	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
24	Endrin	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
25	Endrin aldehyde	

25 Endrin aldehyde

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
25	Endrin aldehyde	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
26	Formaldehyde	Distillation, Colorimetric Method ⁽⁴⁾
27	Free Chlorine	1) Iodometric Method ⁽⁴⁾ 2) DPD Ferrous Titrimetric Method ⁽⁴⁾
28	Heptachlor	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
29	Heptachlor Epoxide	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
30	Hexavalent Chromium	Colorimetric Method ⁽⁴⁾
31	Lead	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
32	Manganese	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
33	Mercury	Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ⁽⁴⁾
34	Methoxychlor	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
35	Nickel	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
36	Oil & Grease	1) Liquid-Liquid, Partition-Gravimetric Method ⁽⁴⁾ 2) Soxhlet Extraction Method ⁽⁴⁾
37	pH	Electrometric Method ⁽⁴⁾
38	Phenols	1) Distillation, Chloroform Extraction Method ⁽⁴⁾ 2) Distillation, Direct Photometric Method ⁽⁴⁾
39	Selenium	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
40	Sulfide	1) Iodometric Method ⁽⁴⁾ 2) Methylene Blue Method ⁽⁴⁾
41	Temperature	Laboratory and Field Methods ⁽⁴⁾
42	Total Dissolved Solids	Dried at 180 °C ⁽⁴⁾
43	Total Kjeldahl Nitrogen	Semi-Micro-Kjeldahl Method ⁽⁴⁾
44	Total Suspended Solids	Dried from 103 to 105 °C ⁽⁴⁾
45	Trivalent Chromium	1) Digestion, Direct Air-Acetylene Flame Method, Colorimetric Method, Calculation ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method, Colorimetric Method, Calculation ⁽⁴⁾
46	Zinc	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾

นำได้ค้น

นำได้ค้น จำนวน 126 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Acenaphthene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
2	Acetone	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
3	Aldrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
4	Anthracene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
5	Antimony	Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
6	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
7	Atrazine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
8	Barium	Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
9	Benz(a)anthracene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
10	Benzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
11	Benzo(b)fluoranthene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
12	Benzo(k)fluoranthene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
13	Benzoic acid	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

14 Benz(a)pyrene

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
14	Benz(a)pyrene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
15	Benz(a,h,i)perylene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
16	Beryllium	Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
17	Bis(2-chloroethyl)ether	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
18	Bis(2-ethylhexyl)phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
19	Bromodichloromethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
20	Bromoform	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
21	Butanol	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
22	Butyl benzyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
23	Cadmium	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ⁽⁴⁾ 3) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
24	Carbazole	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
25	Carbon disulfide	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
26	Carbon tetrachloride	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
27	Chlordane	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
28	p-Chloroaniline	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

29 Chlorobenzene

ลำดับ	สารเคมี	วิธีวิเคราะห์
29	Chlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
30	Chlorodibromomethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
31	Chloroform	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
32	2-Chlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
33	Chromium	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
34	Chromium (III)	1) Digestion, Direct Air-Acetylene Flame Method; Colorimetric Method, Calculation ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method; Colorimetric Method, Calculation ⁽⁴⁾
35	Chromium (VI)	Colorimetric Method ⁽⁴⁾
36	Chrysene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
37	Cyanide	Distillation, Colorimetric Method ⁽⁴⁾
38	2,4-D	Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾
39	DDD	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
40	DDE	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
41	DDT	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
42	Dibenz(a,h)anthracene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

43 Di-n-butyl phthalate

ลำดับ	สารเคมี	วิธีวิเคราะห์
43	Di-n-butyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
44	1,2-Dichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
45	1,3-Dichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
46	1,4-Dichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
47	3,3'-Dichlorobenzidine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
48	1,1-Dichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
49	1,2-Dichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
50	1,1-Dichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
51	cis-1,2-Dichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
52	trans-1,2-Dichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
53	2,4-Dichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
54	1,2-Dichloropropane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
55	1,3-Dichloropropane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
56	1,3-Dichloropropene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
57	Dieldrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
58	Diethyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
59	2,4-Dimethylphenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
60	2,4-Dinitrophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

61 2,4-Dinitrotoluene

ลำดับ	สารเคมี	วิธีวิเคราะห์
61	2,4-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
62	2,6-Dinitrotoluene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
63	Di-n-Octyl phthalate	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
64	Endosulfan	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
65	Endrin	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
66	Ethylbenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
67	Fluoranthene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
68	Fluorene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
69	Heptachlor	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
70	Heptachlor epoxide	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
71	Hexachlorobenzene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
72	Hexachloro 1,3-butadiene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
73	n-Hexane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

74 α -HCH

ลำดับ	สารเคมี	วิธีวิเคราะห์
74	α -HCH	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
75	β -HCH	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
76	γ -HCH	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
77	Hexachlorocyclopentadiene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
78	Hexachloroethane	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
79	Indeno(1,2,3-cd)pyrene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
80	Isophorone	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
81	Lead	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Electrothermal Atomic Absorption Spectrometric Method ⁽⁴⁾ 3) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
82	Manganese	1) Digestion, Direct Air-Acetylene Flame Method ⁽⁴⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽⁴⁾
83	Mercury	Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ⁽⁴⁾
84	Methanol	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
85	Methoxychlor	1) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾
86	Methyl bromide	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽⁴⁾

87 Methylene chloride

ลำดับ	สารเคมี	วิธีวิเคราะห์
87	Methylene chloride	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
88	2-Methylphenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
89	2-Methylnaphthalene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽¹⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
90	Methyl tert-butyl ether	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
91	Naphthalene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽¹⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
92	Nickel	1) Digestion, Direct Air-Acetylene Flame Method ⁽¹⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽¹⁾
93	Nitrobenzene	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
94	N-Nitrosodiphenylamine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
95	N-Nitrosodi-n-propylamine	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
96	Polychlorinated Biphenyls - PCB-1016 - PCB-1221 - PCB-1232 - PCB-1242 - PCB-1248 - PCB-1254 - PCB-1260	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽¹⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
97	Pentachlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
98	pH	Electrometric Method ⁽¹⁾
99	Phenanthrene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽¹⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾

100 Phenol

ลำดับ	สารเคมี	วิธีวิเคราะห์
100	Phenol	1) Distillation, Chloroform Extraction Method ⁽¹⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
101	Pyrene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽¹⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
102	Selenium	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ⁽¹⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽¹⁾
103	Silver	Digestion, Inductively Coupled Plasma Method ⁽¹⁾
104	Styrene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
105	1,1,2,2-Tetrachloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
106	Tetrachloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
107	Toluene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
108	Toxaphene	1) Liquid-Liquid Extraction, Gas Chromatographic Method ⁽¹⁾ 2) Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
109	TPH (C ₅ - C ₆)	1) Purge and Trap, Gas Chromatographic Method ^(1,2,22) 2) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(1,2,21)
110	TPH (C ₁₀ - C ₁₄)	Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(1,2,22)
111	TPH (C ₁₅ - C ₁₉)	Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(1,2,22)
112	1,2,4-Trichlorobenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
113	1,1,1-Trichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
114	1,1,2-Trichloroethane	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
115	Trichloroethylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾

116 2,4,5-Trichlorophenol

ลำดับ	สารเคมี	วิธีวิเคราะห์
116	2,4,5-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
117	2,4,6-Trichlorophenol	Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
118	1,3,5-Trimethylbenzene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
119	Vanadium	Digestion, Inductively Coupled Plasma Method ⁽¹⁾
120	Vinyl acetate	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
121	Vinyl chloride	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
122	m-Xylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
123	o-Xylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
124	p-Xylene	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
125	Xylene (Total)	Purge and Trap Gas Chromatographic/Mass Spectrometric Method ⁽¹⁾
126	Zinc	1) Digestion, Direct Air-Acetylene Flame Method ⁽¹⁾ 2) Digestion, Inductively Coupled Plasma Method ⁽¹⁾

ตารางเคมี (ต่อเนื่องจาก) จำนวน 25 รายการ

ลำดับ	สารเคมี	วิธีวิเคราะห์
1	Antimony	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽¹⁾
2	Arsenic	1) Isokinetic Sampling, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ⁽¹⁾ 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽¹⁾
3	Cadmium	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ⁽¹⁾ 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽¹⁾
4	Carbon Monoxide	Instrumental Analyzer Method ⁽¹⁾
5	Chlorine	Isokinetic Sampling, Ion Chromatographic Method ⁽¹⁾
6	Chromium	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ⁽¹⁾

Chromium (Rp)

ลำดับ	สารเคมี	วิธีวิเคราะห์
6	Chromium (Rp)	2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽¹⁾
7	Cobalt	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽¹⁾
8	Copper	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ⁽¹⁾ 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽¹⁾
9	Cresol	Absorption Sampling, Gas Chromatographic Method ⁽¹⁾
10	Dioxins/Furans	Isokinetic Sampling ⁽¹⁾
11	Hydrogen Chloride	Isokinetic Sampling, Ion Chromatographic Method ⁽¹⁾
12	Hydrogen Fluoride	Isokinetic Sampling, Ion Chromatographic Method ⁽¹⁾
13	Hydrogen Sulfide	Absorption Sampling, Iodometric Method ⁽¹⁾
14	Lead	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ⁽¹⁾ 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽¹⁾
15	Manganese	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ⁽¹⁾ 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽¹⁾
16	Mercury	Isokinetic Sampling, Digestion, Cold Vapor Atomic Absorption Spectrometric Method ⁽¹⁾
17	Nickel	1) Isokinetic Sampling, Digestion, Direct Air-Acetylene Flame Method ⁽¹⁾ 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽¹⁾
18	Opacity	Ringelmann's Method ⁽¹⁾
19	Oxides of Nitrogen	1) Absorption Sampling, Phenoldisulfonic acid Method ⁽¹⁾ 2) Instrumental Analyzer Method ⁽¹⁾
20	Selenium	1) Isokinetic Sampling, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ⁽¹⁾ 2) Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽¹⁾
21	Sulfur Dioxide	1) Absorption Sampling, Barium-Thoron Titrimetric Method ⁽¹⁾ 2) Instrumental Analyzer Method ⁽¹⁾
22	Sulfuric Acid	Isokinetic Sampling, Barium-Thoron Titrimetric Method ⁽¹⁾

23 Total Suspended Particulate

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
23	Total Suspended Particulate	Isokinetic Sampling, Gravimetric Method ⁽³⁾
24	Vanadium	Isokinetic Sampling, Digestion, Inductively Coupled Plasma Method ⁽³⁾
25	Xylene	1) Bag Sampling, Gas Chromatographic Method ⁽³⁾ 2) Adsorption Sampling, Gas Chromatographic Method ⁽³⁾

สิ่งปลูกหรือวัตถุที่พบใช้แล้ว จำนวน 35 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Aldrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
2	Antimony	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
3	Arsenic	1) Waste Extraction, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^(3,6,14) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 3) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^(7,14) 4) Digestion, Inductively Coupled Plasma Method ^(7,14)
4	Barium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
5	Beryllium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
6	Cadmium	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^(3,6,15) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 3) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 4) Digestion, Inductively Coupled Plasma Method ^(7,14)
7	Chlordane	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)

8 Chromium

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
8	Chromium	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^(3,6,15) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 3) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 4) Digestion, Inductively Coupled Plasma Method ^(7,14)
9	Chromium (III)	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method, Waste Extraction, Colorimetric Method, Calculation ^(3,6,13,17) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method, Waste Extraction, Colorimetric Method, Calculation ^(3,6,14,17) 3) Digestion, Flame Atomic Absorption Spectrometric Method, Alkaline Digestion, Colorimetric Method, Calculation ^(7,8,13,17) 4) Digestion, Inductively Coupled Plasma Method, Alkaline Digestion, Colorimetric Method, Calculation ^(7,8,16,17)
10	Chromium (VI)	1) Waste Extraction, Colorimetric Method ^(3,17) 2) Alkaline Digestion, Colorimetric Method ^(8,17)
11	Cobalt	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
12	Copper	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^(3,6,15) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 3) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 4) Digestion, Inductively Coupled Plasma Method ^(7,14)
13	2,4 D	1) Waste Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
14	DDD	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)

15 DDE

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
15	DDE	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
16	DDT	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
17	Dieldrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
18	Endrin	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
19	Heptachlor	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
20	Lead	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^(3,6,15) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 3) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 4) Digestion, Inductively Coupled Plasma Method ^(7,14)
21	Lindane	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
22	Mercury	1) Waste Extraction, Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^(7,19) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 3) Digestion, Cold-Vapor Atomic Absorption Spectrometric Method ^(7,19) 4) Digestion, Inductively Coupled Plasma Method ^(7,14)

Mercury (IIb)

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
22	Mercury (IIb)	5) Thermal Decomposition Amalgamation and Atomic Absorption Spectrometric Method ⁽²⁰⁾
23	Methoxychlor	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
24	Molybdenum	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
25	Nickel	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^(3,6,15) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,6,14) 3) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 4) Digestion, Inductively Coupled Plasma Method ^(7,14)
26	Polychlorinated Biphenyls - Aroclor 1016 - Aroclor 1221 - Aroclor 1232 - Aroclor 1248 - Aroclor 1254 - Aroclor 1260 - 2-Chlorobiphenyl - 2,3-Dichlorobiphenyl - 2,2',5-Trichlorobiphenyl - 2,4',5-Trichlorobiphenyl - 2,2',3,5'-Tetrachlorobiphenyl - 2,2',5,5'-Tetrachlorobiphenyl - 2,3,4,4'-Tetrachlorobiphenyl - 2,2',3,4,5'-Pentachlorobiphenyl - 2,2',4,5,5'-Pentachlorobiphenyl - 2,3,3',4',6-Pentachlorobiphenyl	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(3,9,24) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,24)

Polychlorinated Biphenyls (IIb)

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
27	Polychlorinated Biphenyls(พีบี) - 2,2',3,4,4',5'- Hexachlorobiphenyl - 2,2',3,4,5,5'- Hexachlorobiphenyl - 2,2',3,5,5',6'- Hexachlorobiphenyl - 2,2',4,4',5,5'- Hexachlorobiphenyl - 2,2',3,3',4,4',5'- Heptachlorobiphenyl - 2,2',3,4,4',5,5'- Heptachlorobiphenyl - 2,2',3,4,4',5,6'- Heptachlorobiphenyl - 2,2',3,3',4,4',5,5',6'- Nonachlorobiphenyl Pentachlorophenol	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic/Mass Spectrometric Method ^(3,9,28) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28) Electrometric Method ^(13,32) 3) Waste Extraction, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^(3,9,23) 4) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,8,14) 5) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^(7,21) 6) Digestion, Inductively Coupled Plasma Method ^(7,14) 7) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,8,14) 8) Digestion, Inductively Coupled Plasma Method ^(7,14)
28	pH	
29	Selenium	1) Waste Extraction, Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^(3,9,23) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,8,14) 3) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^(7,21) 4) Digestion, Inductively Coupled Plasma Method ^(7,14)
30	Silver	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,8,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
31	Thallium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,8,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)

32 Toxaphene

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
32	Toxaphene	1) Waste Extraction, Separatory Funnel Liquid-Liquid Extraction, Gas Chromatographic Method ^(7,9,23) 2) Ultrasonic Extraction, Gas Chromatographic Method ^(10,23)
33	Trichloroethylene	1) Waste Extraction, Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(7,12,21) 2) Waste Extraction, Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(7,11,27) 3) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 4) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
34	Vanadium	1) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,8,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
35	Zinc	1) Waste Extraction, Digestion, Flame Atomic Absorption Spectrometric Method ^(3,8,15) 2) Waste Extraction, Digestion, Inductively Coupled Plasma Method ^(3,8,14) 3) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 4) Digestion, Inductively Coupled Plasma Method ^(7,14)

สิ้น จำนวน 125 รายการ

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
1	Acenaphthene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
2	Acetone	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
3	Aldrin	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
4	Anthracene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25)

Anthracene (พีบี)

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
4	Anthracene (พีบี)	2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
5	Antimony	Digestion, Inductively Coupled Plasma Method ^(7,14)
6	Arsenic	1) Digestion, Hydride Generation/Atomic Absorption Spectrometric Method ^(7,14) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
7	Atrazine	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
8	Barium	Digestion, Inductively Coupled Plasma Method ^(7,14)
9	Benz(a)anthracene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
10	Benzene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
11	Benzo(b)fluoranthene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
12	Benzo(k)fluoranthene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
13	Benzoic acid	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
14	Benzo(a)pyrene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
15	Benzo(g,h,i)perylene	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
16	Beryllium	Digestion, Inductively Coupled Plasma Method ^(7,14)

17 Bis(2-chloroethyl)ether...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
17	Bis(2-chloroethyl)ether	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
18	Bis(2-ethylhexyl)phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
19	Bromodichloromethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
20	Bromoform	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
21	Butanol	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
22	Butyl benzyl phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
23	Cadmium	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 2) Digestion, Inductively Coupled Plasma Method ^(7,14)
24	Carbazole	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
25	Carbon disulfide	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
26	Carbon tetrachloride	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
27	Chlordane	1) Ultrasonic Extraction, Gas Chromatographic Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
28	p-Chloroaniline	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
29	Chlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
30	Chlorodibromomethane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
31	Chloroform	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(13,27)
32	2-Chlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)

33 Chromium...

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
33	Chromium	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(1,15) 2) Digestion, Inductively Coupled Plasma Method ^(9,14)
34	Chromium (III)	1) Digestion, Flame Atomic Absorption Spectrometric Method, Alkaline Digestion, Colorimetric Method, Calculation ^(7,8,15,17) 2) Digestion, Inductively Coupled Plasma Method, Alkaline Digestion, Colorimetric Method, Calculation ^(7,8,14,17)
35	Chromium (VI)	Alkaline Digestion, Colorimetric Method ^(8,17)
36	Chrysene	1) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
37	Cyanide	Extraction, Distillation, Colorimetric Method ^(29,30)
38	2,4-D	Ultrasonic Extraction, Gas Chromatographic Method ⁽²⁶⁾
39	DDD	1) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
40	DDE	1) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
41	DDT	1) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
42	Dibenz(a,h)anthracene	1) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
43	Di-n-butyl phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
44	1,2-Dichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(12,27)

45 1,3-Dichlorobenzene

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
45	1,3-Dichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(12,27)
46	1,4-Dichlorobenzene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(12,27)
47	3,3'-Dichlorobenzidine	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
48	1,1-Dichloroethane	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(12,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
49	1,2-Dichloroethane	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(12,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
50	1,1'-Dichloroethylene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(12,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
51	cis-1,2-Dichloroethylene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(12,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
52	trans-1,2-Dichloroethylene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(12,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
53	2,4-Dichlorophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
54	1,2-Dichloropropane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(12,27)
55	1,3-Dichloropropane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(12,27)
56	1,3-Dichloropropene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(12,27)
57	Dieldrin	1) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)

58 Diethyl phthalate

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
58	Diethyl phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
59	2,4-Dimethylphenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
60	2,4-Dinitrophenol	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
61	2,4-Dinitrotoluene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
62	2,6-Dinitrotoluene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
63	Di-n-Octyl phthalate	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
64	Endosulfan	1) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
65	Endrin	1) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
66	Ethylbenzene	1) Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(12,27) 2) Equilibrium Headspace, Gas Chromatographic/Mass Spectrometric Method ^(11,27)
67	Fluoranthene	1) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
68	Fluorene	1) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
69	Heptachlor	1) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
70	Heptachlor epoxide	1) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,23)

Heptachlor epoxide (คป)

ลำดับ	สารมลพิษ	วิธีวิเคราะห์
70	Heptachlor epoxide (คป)	2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
71	Hexachlorobenzene	1) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
72	Hexachloro-1,3-butadiene	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(12,27)
73	n-Hexane	Purge and Trap, Gas Chromatographic/Mass Spectrometric Method ^(12,27)
74	α-HCH	1) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
75	β-HCH	1) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
76	γ-HCH	1) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,23) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
77	Hexachlorocyclopentadiene	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
78	Hexachloroethane	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
79	Indeno(1,2,3-cd)pyrene	1) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,25) 2) Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
80	Isophorone	Ultrasonic Extraction, Gas Chromatographic/Mass Spectrometric Method ^(10,28)
81	Lead	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 2) Digestion, Inductively Coupled Plasma Method ^(9,14)
82	Manganese	1) Digestion, Flame Atomic Absorption Spectrometric Method ^(7,15) 2) Digestion, Inductively Coupled Plasma Method ^(9,14)

83 Mercury...

ເອກສານອ້າງອີງ

14 United States

14. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Inductively Coupled Plasma-Optical Emission Spectrometry. SW-846 Method 6010D, 2014.

15. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Flame Atomic Absorption Spectrophotometry. SW-846 Method 7000B, 2007.

16. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Arsenic (Atomic Absorption, Gaseous Hydride). SW 846 Method 7061A, 1992.

17. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Chromium, Hexavalent (Colorimetric). SW-846 Method 7196A, 1992.

18. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Mercury in Liquid Waste (Manual Cold Vapor Technique). SW-846 Method 7470A, 1994.

19. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Mercury in Solid or Semisolid Waste (Manual Cold-Vapor Technique). SW-846 Method 7471B, 1998.

20. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Mercury in Solids and Solutions by Thermal Decomposition, Amalgamation, and Atomic Absorption Spectrophotometry. SW-846 Method 7473, 2007.

21. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Selenium (Atomic Absorption, Borohydride Reduction). SW-846 Method 7742, 1994.

22. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Nonhalogenated Organics Using GC/FID. SW-846 Method 8015D, 2003.

23. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Organochlorine Pesticides by Gas Chromatography. SW-846 Method 8081B, 2007.

24. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Polychlorinated Biphenyls (PCBs) by Gas Chromatography. SW-846 Method 8082A, 2007.

25. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Organochlorine Pesticides by Gas Chromatography. SW-846 Method 8081B, 2007.

26. United States Environmental Protection Agency. Test Methods for Evaluation Solid Waste Physical/Chemical Methods. Chlorinated Herbicides by GC Using Methylation or Pentafluorobenzoylation Derivatization. SW-846 Method 8151A, 1996.

27. United States

9.52

สำเนาถูกต้อง R/h



UNITED ANALYST AND ENGINEERING
CONSULTANT COMPANY LIMITED

บริษัท ยูไนเต็ด แอนนาลิสต์ แอนด์ เอ็นจิเนียริง คอนซัลแตนท์ จำกัด
เลขที่ 3 ซอยอุดมสุข 41 ถนนสุขุมวิท แขวงบางจาก เขตพระโขนง กรุงเทพฯ 10260
โทรศัพท์ 0-2763-2828 E-mail address : uae@uaeconsultant.com